



Discover Dexcom G5 Mobile Continuous Glucose Monitoring (CGM) System


CGM unveils what no meter can. While meters give a number for a single point in time, CGM provides dynamic glucose information by showing where your glucose is, where it's going, and how fast it's getting there. Better yet, with Dexcom G5 Mobile, you can eliminate the need for confirmatory fingersticks, for daily glucose management decisions.*

*If your glucose alerts and readings do not match your symptoms or expectations or you are taking medications containing paracetamol/acetaminophen, you should perform a fingerstick to confirm your

blood glucose level. A fingerstick is needed every 12 hours to calibrate.

Contact Your Dexcom Representative >



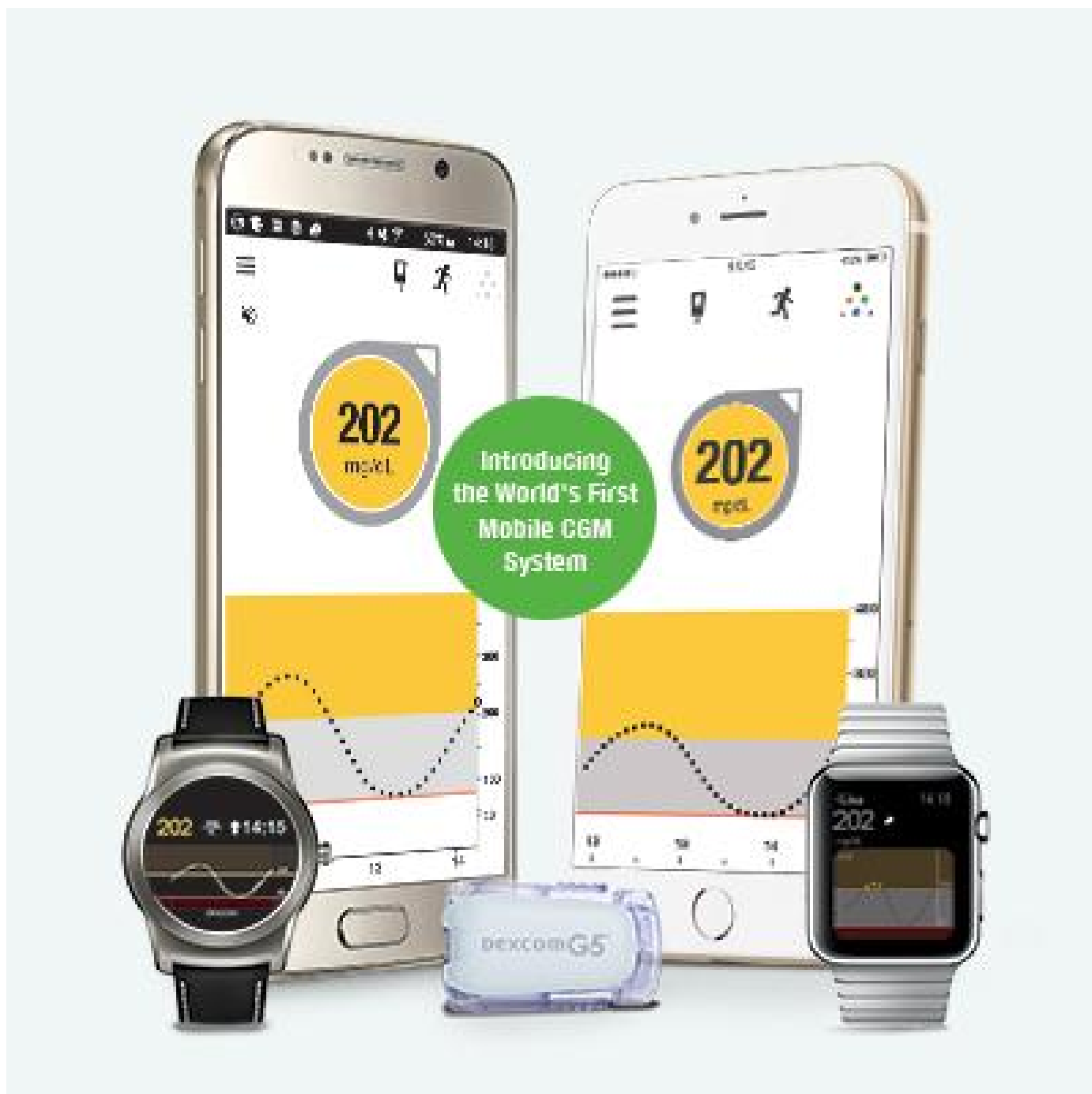

Meet the Dexcom G5 Mobile CGM System



1. Small Sensor that measures glucose levels just underneath the skin.
2. Transmitter that is fastened on top of the sensor and sends data wirelessly to your compatible smart device or your receiver.
3. A Compatible Smart Device with the Dexcom G5 Mobile app. View your glucose trends in vivid colors to know when it's high, low or within range.

A Dexcom G5 Mobile Receiver can also be acquired as an optional display device.






Compatible with Android and Apple Smart Devices³

Monitor your glucose on the go with the world's first mobile CGM platform

Dexcom G5 Mobile is the world's first real time continuous glucose monitoring (CGM) system approved for adults and children 2 years and older. Using Bluetooth technology that is built into the transmitter, Dexcom G5 Mobile allows for complete remote viewing of glucose levels, trends and data from your [compatible smart device](#).



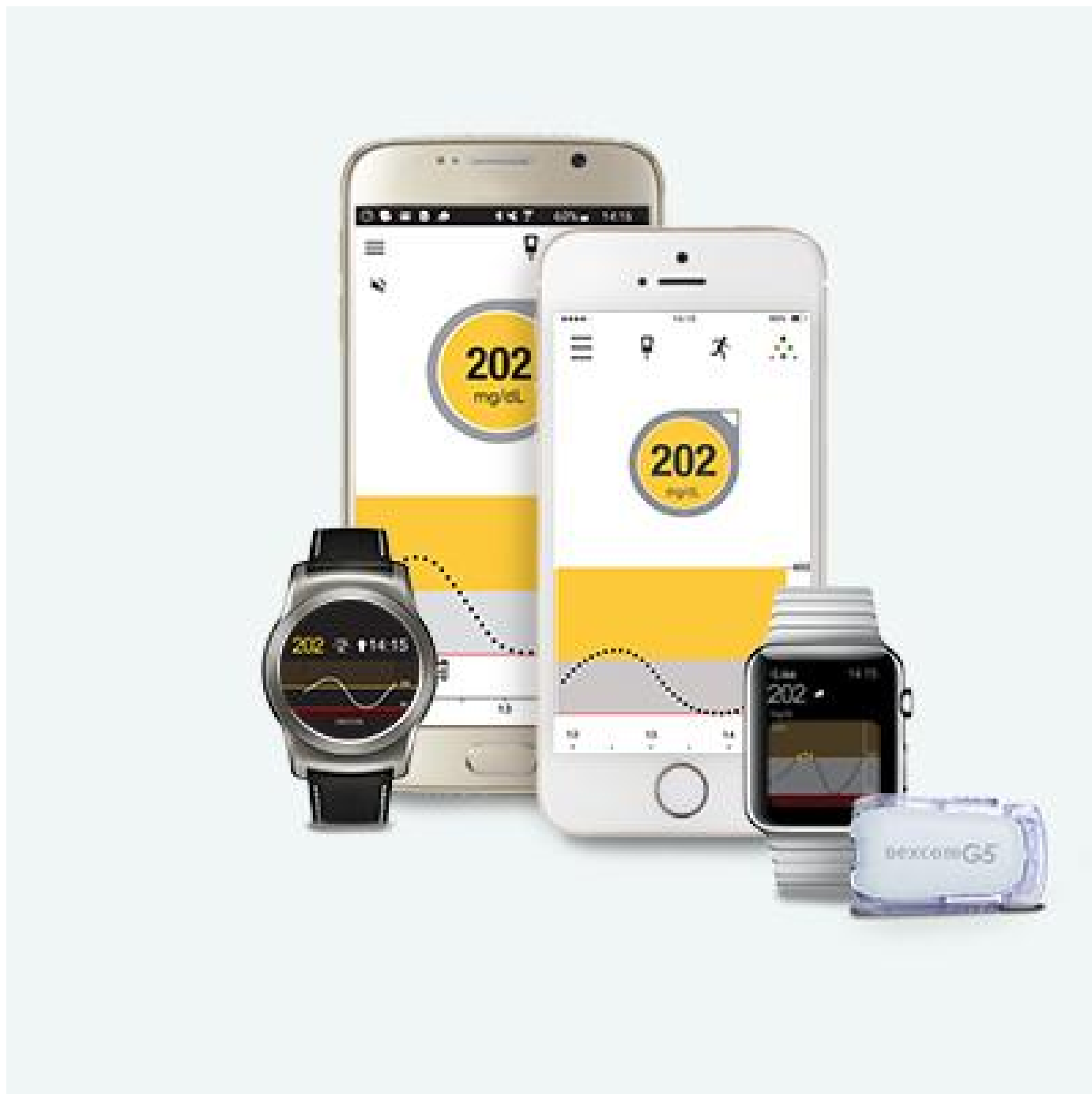

Receive High and Low Glucose Alerts



With the Dexcom G5 Mobile CGM System, you'll be alerted directly on your compatible smart device when you're heading too high or too low. Alerts and alarm sounds can also be customized to appear as a text message, allowing for additional discretion and privacy.






Track Your Actions

Now you can keep track of events that affect your glucose levels on your compatible smart device. The Dexcom G5 Mobile CGM

System provides a platform to enter customizable events, giving you the ability to track how your daily activities influence your glucose trends.

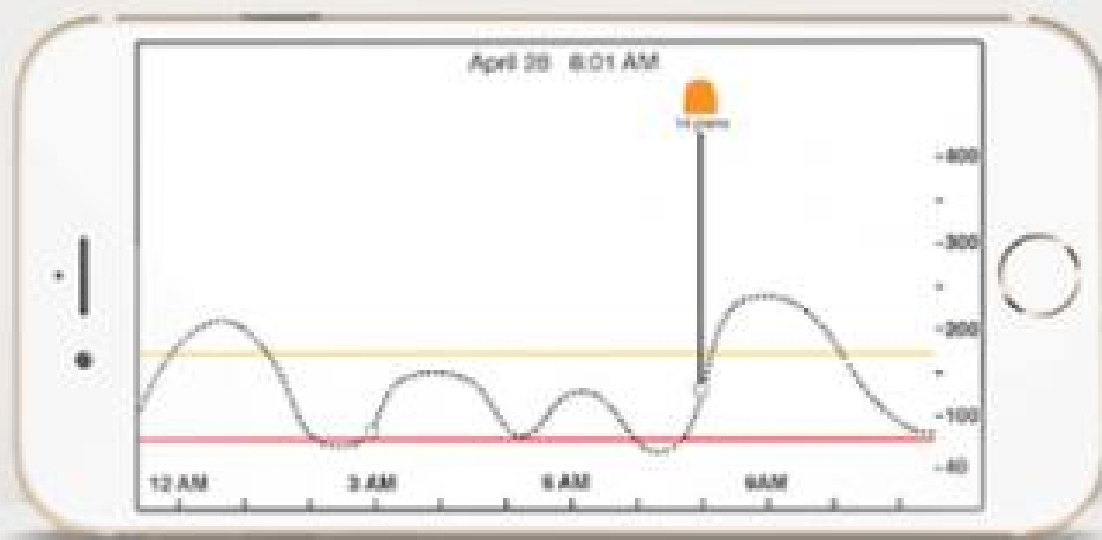

Share your CGM data with Dexcom Share

With the Dexcom G5 Mobile CGM System, you can stay connected in real-time to those nearest and dearest to you - all from your [compatible smart device](#).

To get started, just follow these 3 simple steps:

1. Enable the Dexcom Share feature from your G5 Mobile App
2. Invite up to 5 people to Follow you
3. Followers will only need to download the Dexcom Follow App and have an internet connection

With real-time glucose readings appearing directly on your Followers smart device,

you'll be able to share your glucose levels with your circle of support to help you better manage your diabetes and have peace of mind.





Accuracy and Sensing Performance

When it comes to glucose sensing, Dexcom sets the standard. Our outstanding accuracy, along with exceptional sensing performance, results in a reliable CGM you can trust. With the Dexcom G5 Mobile CGM System, you'll get a CGM that's not only appealing but one that delivers results consistently and accurately^{1,2}.

When **LOWER** is better

Adults ^{1*}	Performance Metrics	Pediatrics ^{2*}
9%	Overall Accuracy Mean ARD% (MARD), 40-400 mg/dL (2.22-22.22 mmol/L) (% average absolute error versus reference across all glucose levels)	10%
7%	Sensor to Sensor Variation Percent Coefficient of Variation, CV	7%

Adults ^{1*}	Performance Metrics	Pediatrics ^{2*}
	(Consistent performance between each sensor)	
Day 1: 11%	Accuracy Over Time Mean ARD% (MARD), 40-400 mg/dL (2.22-22.22 mmol/L) (Consistent sensor performance over time)	Day 1: 13%
Day 4: 8%	Average Daily Calibration 2X/day	Day 4: 8%
Day 7: 9%		Day 7: 10%

When **HIGHER** is better

Adults ^{1*}	Performance Metrics	Pediatrics ^{2*}
92% [99.5%]	Clinical Accuracy % of readings that were in the Clarke Error Grid (CEG) A Zone [% CEG A+B Zone]	90% [98.5%]

Adults ^{1*}	Performance Metrics	Pediatrics ^{2*}
98%	Sensor Life % of sensors that last up to 7 days use	94%

* Reference is YSI (Yellow Springs Laboratory Instrument)

References :

1. Bailey, Chang, Christiansen Clinical Accuracy of a Continuous Glucose Monitoring System with an Advanced Algorithm J Diabetes Sc Tech 2015 Vol 9(2) 209-214
2. Laffel Improved Accuracy of Continuous Glucose Monitoring Systems in Pediatric Patients with Diabetes Mellitus - Results from Two Studies Diabetes Technology & Therapeutics Volume 18, Supplement 2, 2016
3. [Click here](#) for a list of compatible smart devices.

Contact Information

Gulf Pharmaceutical Industries (Julphar)
Scientific Office
Building 25



Dubai Healthcare City
Tel: +971 4 4385 222
Direct: +971 56 509 4803
Fax: +971 4 4385 252
dexcom@julphar.net
www.julphardiabetes.net

Downloads & Guides (Arabic)

Dexcom G5 Mobile CGM System (Arabic)

[Dexcom G5 Mobile User's Guide](#)

[Dexcom G5 Mobile Quick Start Guide](#)

[Dexcom G5 Mobile Tutorial](#) (English Only)

[Dexcom Share/Follow User Guide](#)

Dexcom G4 PLATINUM (Arabic)

[Dexcom G4 Platinum User's Guide mmol/L](#)

[Dexcom G4 Platinum User's Guide mg/dL](#)

[Dexcom G4 Platinum Quick Start Guide mmol/L](#)

[Dexcom G4 Platinum Quick Start Guide mg/dL](#)

[Dexcom G4 Platinum Tutorial](#)



Downloads & Guides (English)

Dexcom G5 Mobile CGM System (English)

[Dexcom G5 Mobile User's Guide](#)

[Dexcom G5 Mobile Quick Start Guide](#)

[Dexcom G5 Mobile Tutorial](#)

[Dexcom Share/Follow User Guide](#)

Dexcom G4 PLATINUM (English)

[Dexcom G4 Platinum User's Guide mmol/L](#)

[Dexcom G4 Platinum User's Guide mg/dL](#)

[Dexcom G4 Platinum Quick Start Guide mmol/L](#)

[Dexcom G4 Platinum Quick Start Guide mg/dL](#)

[Dexcom G4 Platinum Tutorial](#)

Dexcom STUDIO (English Only)

[Dexcom Studio Software Download](#)

[Dexcom Studio User's Guide](#)



Dexcom Studio Quick Start Guide

Terms of Use
Privacy Policy
Trust Center

Safety
Information

Dexcom, Dexcom Clarity, Dexcom Follow, Dexcom One, Dexcom Share, Share are trademark or registered trademarks in the U.S. and may be in other countries. © 2023 Dexcom, Inc. All rights reserved.
LBL015482 Rev001

