



User Guide

Guardian[®] REAL-Time

CONTINUOUS GLUCOSE MONITORING SYSTEM



DexCom Exhibit 2008
Abbott v. DexCom
IPR2024-00840

Guardian[®] REAL-Time
Continuous Glucose Monitoring System
User Guide

©2006, Medtronic MiniMed, Inc. All rights reserved.

This product is covered by U.S. Patent No. 6,809,653. Other U.S. and/or foreign patents may be pending.

Guardian®, Medtronic CareLink® and Sen-serter® are registered trademarks of Medtronic MiniMed, Inc.

CareLink Personal™, ComLink™, and Medtronic MiniLink™ are trademarks of Medtronic MiniMed, Inc.

Energizer® is a registered trademark of Eveready Battery Company.

Glucagon Emergency Kit® is a registered trademark of Eli Lilly and Company.

Cavicide® is a registered trademark of Metrex Research Corporation, an indirect subsidiary of Sybron Dental Specialties, Inc.

LifeScan®, OneTouch®, OneTouch® Ultra® and OneTouch® UltraLink™ are trademarks of LifeScan, Inc.



Medtronic

Medtronic MiniMed
Northridge, CA 91325
USA
800-646-4633
818-576-5555
www.medtronicdiabetes.com

R_x Only



6025254-015 070408

REF CSS72, CSS7100, CSS7100K, CSS7707

Contacts:**Africa:**

Medtronic Africa (Pty) Ltd.
Tel: +27 (0) 11 677 4800

Australia:

Medtronic Australasia Pty. Ltd.
Tel: 1800 668 670 (product orders)
Tel: 1800 777 808 (customer help)

Belgie/Belgique:

N.V. Medtronic Belgium S.A.
Tel: 02-456-0900

Brasil:

Medtronic Comercial Ltda.
Tel: +(11) 3707-3707

Canada:

Medtronic of Canada Ltd.
Tel: 1-800-284-4416 (toll free/sans-frais)

Ceska Republika:

Medtronic Czechia s.r.o.
Tel: +420 296 579 580

China:

Medtronic (Shanghai) Ltd.
Tel: +86 40 0820 1981
or 80 0820 1981

Danmark:

Medtronic Danmark A/S
Tel: +45 32 48 18 00

Deutschland:

Medtronic GmbH
Geschäftsbereich Diabetes
Telefon: +49 2159 8149-370
Telefax: +49 2159 8149-110
24-Std-Hotline: 0800 6464633

East Asia:

Medtronic International Ltd.
Tel: +852 2891 4300

Eire:

Accu-Science LTD.
Tel: +353 45 433000

España:

Medtronic Ibérica S.A.
Tel: + 34 91 625 05 42
Fax: +34 91 625 03 90
24 hour: +34 901 120 335

Europe:

Medtronic Europe S.A. Europe, Middle East and Africa Headquarters
Tel: +41 (0) 21-802-7000

France:

Medtronic France S.A.S.
Tel: +33 (0) 1 55 38 17 00

Hellas:

Medtronic Hellas S.A.
Tel: +30 210677-9099

Italia:

Medtronic Italia S.p.A.
Tel. +390 2 24137 261
Servizio assistenza tecnica: No verde
24h: 800 712 712
Tel. 06-328141
Fax 06-3215812

Japan:

Medtronic Japan Co. Ltd.
Tel: +81-3-6430-2019

Latin America:

Medtronic, Inc.
Tel. (1305)-500-9328
Fax (1786)-709-4244

Latvija:

Ravemma Ltd.
Tel: +371 7273780

Magyarország:

Medtronic Hungária Kft.
Tel: +36 1 889 0600

Malaysia:

Medtronic International Ltd.
Tel: (+65) 6436-5094
or (+65) 6436-5097

Middle East and North Africa:

Regional Office
Tel: +961-1-370 670

Nederland, Luxembourg:

Medtronic B.V.
Tel: +31 (0) 45-566-8290

New Zealand:

Medica Pacifica
Tel: +0800 106 100
24-hr After-Hours:
+0800 633 487

Norge:

Medtronic Norge A/S
Tel: +47 67 10 32 00
Fax: +47 67 10 32 10

РОССИЯ:

Medtronic B. V.
Tel: +7 095 5807377/1161

Polska:

Medtronic Poland Sp. z.o.o.
Tel: +48 22 4656 900

Portugal:

Medtronic Portugal Lda
Tel: +351 21 7245100
Fax: +351 21 7245199

Puerto Rico:

Medtronic Puerto Rico
Tel: 787-753-5270

Republic of Korea:

Medtronic Korea, Co., Ltd.
Tel: +82.2.3404.3600

Schweiz:

Medtronic (Schweiz) AG
Tel: + 41 (0)31 868 0160
24-Stunden-Hotline 0800 633 333
Fax Allgemein: +41 (0)31 8680199

Serbia & Montenegro:

Novolab Serbia & Montenegro
Epsilon Research Ltd.
Tel: +381 63 219 827

Singapore:

Medtronic International Ltd.
Tel: (+65) 6436-5094
or (+65) 6436-5097

Slovenija:

Zaloker & Zaloker d.o.o.
Tel: +34 91 625 0400

Slovenská Republika:

Medtronic Slovakia o.z.
Tel: +421 2 68 20 69 12

Suomi:

Medtronic Finland OY
Tel: +358-9-755 25 00

Sverige:

Medtronic AB
Tel: +46 8 568 585 00
Fax: +46 8 568 585 01

Taiwan:

Medtronic-Mediland Ltd.
Tel: +886.2.2517.0308

Thailand:

Medtronic (Thailand) Ltd.
Tel: (+65) 6436-5094
or (+65) 6436-5097

Türkiye:

Medtronic Medikal Teknoloji
Ticaret Ltd. Şirketi.
Tel: +90 216 4694330

USA:

Medtronic Diabetes Global
Headquarters
Tel: 800-826-2099
24-hour Help Line: 818-576-5555
To order supplies 800-843-6687

United Kingdom:

Medtronic Ltd.
Tel: +44 1923-205167

Yisrael:

Agentek
Tel: +972 3649 3111

Österreich:

Medtronic Österreich GmbH
Tel: +43 (0) 1 240 44-0
24 - Stunden - Hotline 0820 820 190

Contents

Chapter 1: Introduction	1
Assistance	2
Emergency kit	2
Your system	3
Accessories	4
Conventions used in this guide	5
User safety	6
Indications: CSS7100	6
Indications: CSS7100K (pediatric version)	6
Contraindications	7
Safety Warnings	7
Sensor	7
Transmitter	7
Meter	7
Precautions	7
Avoid extreme temperatures	7
Sensor	8
Meter	8
Adverse reactions	8
Damage due to radiation	8
Compliance information	9
RF interference from other devices	10
Chapter 2: The basics	11
The monitor	11
Installing the battery	12
The monitor buttons	15

An important message about alarms and alerts	16
The monitor screen	17
The HOME screen	17
Icons	17
Sensor icon	17
Battery icon	18
Time of day icon	18
Alert icon	18
Alarm icon	18
Scroll bar	18
Screen backlight	18
Alarm/alert types	19
Operating modes	19
Normal mode	19
Special mode (Alerts)	19
Attention mode (Alarms)	20
Menus	20
MAIN MENU	20
SENSOR MENU	20
CAPTURE EVENT	20
UTILITIES MENU	21
The STATUS screen	22

Chapter 3: Entering your settings 23

Basic steps	23
Setting the time and date	24
Selecting the language	26
Setting your Alarm/Alert Type	26
Entering your glucose monitoring settings	27
Turning on the sensor	27
Selecting the BG units	28
Turning on the Glucose Alerts	28
Setting the Glucose Limits	28
Setting the Glucose Limits start time	30
Setting the High Snooze	31
Setting the Low Snooze	32

Setting the Predictive glucose alert	32
Setting the Rate Alerts	34
Setting the AUC Limits	35
Setting the Alarm Snooze	35
Setting the Cal Reminder	36
Entering the transmitter identification number	36
Setting up the Missed Data option	37
Setting the sensor glucose graph timeout	37
Reviewing your settings	38
Using the Sensor Demo	38

Chapter 4: Setting up your transmitter, sensor and meter 39

Basic steps	39
Setting up the transmitter	40
MiniLink charger	40
Installing a new charger battery	40
Charging the transmitter	41
Inserting the sensor	42
Before you begin	43
Make sure that the transmitter is fully charged	43
Make sure that the sensor is ready to be inserted	43
Where to insert the sensor	43
Do NOT insert the sensor in	43
How to prepare the insertion site	44
Connecting the transmitter to the sensor	46
Starting the sensor	47
Entering the first meter blood glucose	47
Meter option	48
Meter rules	49
Add, delete, review meter IDs	50

Chapter 5: Using your system 51

Calibrating your system	51
Guidelines for calibrating the system	52
Capturing Events	54
Entering non-calibration meter BG measurements	55

Entering insulin injection information	55
Entering information about carbohydrates	56
Entering exercise information	57
Entering Other markers	58
Reading the sensor glucose graphs	58
Opening and viewing the graphs	59
The graphs	60
Examples of real-time sensor glucose graphs	61
3-hour graph	61
6-hour graph	62
12-hour graph	63
24-hour graph	64
How to check for rapid changes in sensor glucose	65
The rapid change arrows	65
Checking the status of the system	67
Checking the monitor battery status	67
Checking the sensor status	68
Checking the transmitter status	69
How to make glucose alerts silent	70
Viewing the sensor update history	71
Viewing your current settings	72
Disconnecting the transmitter and removing the sensor	72
Disconnecting the sensor from the transmitter	72
Removing the sensor	73
Using your system in water	73
Chapter 6: Utilities	75
The data management feature	75
Sensor glucose information	76
Meter blood glucose information	77
The Alarm Clock	78
Reviewing your current Alarm Clock reminders	79
Deleting Alarm Clock reminders	79
The Block feature	80
The Lock Keypad feature	80
Running the Selftest	81

User settings	82
Saving your settings	82
Restoring your settings	83
Clearing your settings	83
Viewing your User Settings history	84
Chapter 7: Therapy software	85
CareLink Personal™	85
Reports	85
Logbook	86
Chapter 8: Troubleshooting and alarms/alerts	87
About alarms and alerts	87
What to do when you get an alarm/alert	88
Sensor alerts	89
WEAK SIGNAL	89
LOST SENSOR	89
LOW TRANSMTR	89
CHARGE TRANSMTR	90
CHANGE SENSOR	90
SENSOR END	90
CAL ERROR	90
METER BG NOW	91
METER BG BY	91
LOW SG	91
HIGH SG	91
LOW PREDICTED	91
HIGH PREDICTED	91
FALL RATE	92
RISE RATE	92
SENSOR ERROR	92
LOW BATTERY	92
Viewing your sensor alert history	93
System alarms	94
A (Alarm) (nn)	94
BATT OUT LIMIT	94

BUTTON ERROR	94
CHECK SETTINGS	94
E (Error) (nn)	94
FAILED BATT TEST	95
OFF NO POWER	95
RESET	95
WEAK BATTERY	95
Viewing your system alarm history	95
Troubleshooting the MiniLink transmitter and charger	97
Checking battery status of the charger and transmitter	97
Troubleshooting the transmitter	99
Connecting the tester	99
Disconnecting the tester	100
Troubleshooting the monitor	100
My monitor will not display my BG measurement from my meter	100
Reconnect old sensor	101
Find lost sensor	101
What happens if I leave the monitor battery out for more than ten minutes?	101
Why doesn't my monitor battery last very long?	102
What is a CHECK SETTINGS alarm?	102
My screen looks distorted	102
I dropped my monitor	102
I submerged my monitor in water	103
I cannot get to the User Settings screen	103

Chapter 9: System maintenance 105

Replacing the battery	105
Storing the system	105
Storing the monitor	105
Storing the MiniLink transmitter	105
Storing the Sen-serter	105
Cleaning the system	106
Cleaning the monitor	106
Cleaning the Sen-serter	107
Cleaning the transmitter	107
Cleaning the transmitter	108

Chapter 10: System specifications	109
Alarm/alert types	110
Alarm/alert messages and codes	110
Alarm/alert history	110
Backlight	110
Screen and menu timeout	110
Sensor glucose graph timeout	111
High and low Glucose Limits (CSS7100)	111
High and low Glucose Limits (CSS7100K)	111
Daily totals	112
Default screen	112
Default settings	112
ENTER BG meter values	119
Meter ID entries	119
Power supply	119
System safety checks	119
Monitor size	119
Monitor weight	120
TIME/DATE SET screen	120
STATUS screen	120
Guidance and manufacturer's declaration	121
MiniLink specifications	127
Symbol table	128
Appendix A: Sensor accuracy	131
Performance results in adults	132
Accuracy of Guardian RT readings	133
Precision of Guardian RT readings	137
Low and High Alerts in adults	138
The Low Glucose Alert	138
The High Glucose Alert	139
Guardian RT Sensor Performance and Calibration Stability As a Function of Time	141
Effects of calibration frequency	144
Time-Elapsed Plots	145
Performance results in children and adolescents	148
Accuracy of Guardian RT readings	148

Low and High alerts in children and adolescents	152
The Low Glucose Alert	152
The High Glucose Alert	153
Glossary	155
Index	161

Introduction

Thank you for choosing Medtronic Diabetes (legally known as Medtronic MiniMed, Inc.) as your diabetes management partner.

Your new Guardian® REAL-Time Continuous Glucose Monitoring (CGM) System is a big step forward in diabetes management for people with diabetes mellitus because it:

- provides continuous or periodic monitoring of glucose levels in the fluid under your skin;
- provides real-time glucose values so that you can track glucose concentration patterns and possibly identify episodes of low and high blood glucose;
- stores glucose data so that it can be analyzed to track patterns and/or downloaded to PC software for analysis of historical glucose values; and
- alerts you in real time if a glucose level falls below or rises above pre-set values.

Pre-set values are not intended to be used directly for making therapy adjustments, but rather to provide an indication of when a meter blood glucose measurement may be required. All therapy adjustments should be based on measurements obtained using a home glucose meter and not on Guardian REAL-Time CGM System values. Therefore, we strongly recommend that you work closely with your healthcare professional when you start your Guardian REAL-Time CGM System, and for ongoing support.

In this chapter:

- Assistance
- Emergency kit
- Your system
- Accessories
- Conventions used in this guide
- User Safety
- Compliance information

Assistance

Medtronic Diabetes provides a 24 Hour HelpLine for assistance. The HelpLine is staffed with representatives who are trained in the setup and operation of the Guardian REAL-Time CGM System and are able to answer related questions. When calling the HelpLine, please have your Guardian monitor available. The monitor serial number and the 24 Hour HelpLine phone number are on the back of your monitor.

Department	Telephone number
24 Hour HelpLine (calls within the United States)	800.646.4633
24 Hour HelpLine (calls outside the United States)	+1.818.576.5555
Web site	www.medtronicdiabetes.com

Emergency kit

You should keep a diabetes emergency kit with you at all times - and be sure a family member, friend or co-worker knows where it is. Your emergency kit may include these items:

- Fast-acting glucose tablets and Glucagon Emergency Kit[®]
- Urine ketone monitoring supplies and blood glucose monitoring supplies
- Insulin syringe and insulin (with dosage instructions from your healthcare professional)
- Guardian REAL-Time CGM System Wallet Card and extra AAA alkaline batteries for your monitor (Energizer[®] brand is recommended)
- Transparent dressing

2 Introduction

Your system

This section describes the parts of your Guardian REAL-Time CGM System.

- **Monitor** — Your monitor (CSS7100/CSS7100K) shows real-time glucose measurements, tells you about changes in glucose levels, provides historical information, and displays system information such as low battery alerts.



- **Transmitter** — The Medtronic MiniLink™ transmitter (MMT-7703) is a small rechargeable device that connects to the glucose sensor that is inserted into the body. When connected to the sensor, the MiniLink transmitter automatically initializes the sensor, collects glucose data, and wirelessly sends the data to your monitor. It comes with a blue tester (MMT-7706) and a charger (MMT-7705).



Transmitter



Tester



Charger

- **Sensor** — The sensor (MMT-7002/MMT-7003) is a membrane-covered electrode inserted through the skin into the fatty layer underneath with a device called the Sen-serter®. The sensor tells the transmitter how much glucose is in the interstitial fluid at the insertion site. The transmitter sends this data to the



monitor and the monitor screen shows you the reading.

- **CareLink USB Upload Device** — The Medtronic Diabetes CareLink USB (MMT-7305) is a radio frequency (RF) device that plugs into your computer's USB port and uploads Guardian REAL-Time CGM System data to your diabetes management software.

Accessories

Medtronic offers the following accessories to help you better use and more easily wear your monitor.

- **Meter** — Your Guardian REAL-Time can be used with an optional blood glucose meter powered by MWT1 technology (where or if available) that you can program to automatically transmit blood glucose test results to your Guardian monitor. MWT1 is the wireless Radio Frequency (RF) technology that is used to transmit information from the meter to the monitor. (You can also use a non-RF meter with your monitor, and manually enter your BG measurements.)
- **Holster** — To wear the monitor on your belt.
- **Monitor clip** — To wear the monitor underneath your clothing.
- **Leather case** — Fine leather lined with nylon. Styling complements business and formal wear. Velcro flap provides easy programming access. Wear it vertically with the built-in belt clip.

To order supplies, call 800-646-4633, +1-828-362-5958 (outside U.S.), refer to the contacts list at the beginning of this user guide, or visit our web site at www.medtronicdiabetes.com.

Conventions used in this guide

The conventions used in this guide are defined in the following table.

Conventions	Meaning
Press	Push and release the button.
Hold	Press and keep pressure on the button.
Press ESC, ACT	Press the ESC button, then press the ACT button. When two buttons are separated by a comma (,) it means to press the first button, then press the next button.
Select	Choose an item or a value from a menu or screen. Use the  and  buttons to make your selections.
Button names	Always bold and uppercase; for example, UP , DOWN , ACT .
Screen and menu names	Always uppercase; for example, MAIN MENU, STATUS screen.
Menu items	Always bold; for example, 12-Hour Setup , On , Off
CAUTION	A potential hazard that, if not avoided, may cause damage to your Guardian REAL-Time CGM System.
WARNING	A potential hazard that, if not avoided, may cause death or serious injury. It may also describe potential serious adverse reactions and safety hazards.
Menu item selections	The items on a menu you need to select to go to another menu or a screen. Each menu item is separated by the (>) symbol. For example: MAIN MENU > Utilities > Alarm The actual steps used to complete menu item selection shown above are: 1 Open the MAIN MENU. 2 Select Utilities, press ACT. 3 Select Alarm, press ACT. To make the steps for getting to a specific menu or screen easier to understand, the menu item selections do not include the steps for pressing the ACT button.

User safety

This section includes important safety information such as indications, contraindications, safety warnings, potential adverse reactions, and how to protect the system from radiation exposure damage.

Indications: CSS7100

The Guardian REAL-Time CGM System (CSS7100) is indicated for continuous or periodic monitoring of glucose levels in the fluid under the skin, in adults (18 years and older) with diabetes mellitus, for the purpose of improving diabetes management. It alerts if a glucose level falls below, or rises above, pre-set values. Values are not intended to be used directly for making therapy adjustments, but rather to provide an indication of when a meter blood glucose measurement may be required. All therapy adjustments should be based on measurements obtained using a home glucose meter and not on Guardian REAL-Time CGM System values.

The Guardian REAL-Time CGM System provides real-time glucose values that allow users to track patterns in glucose concentrations and to possibly identify episodes of low and high blood glucose. It also stores the data so that it can be analyzed to track patterns. Glucose data can be further downloaded to PC software for analysis of historical glucose values. The Guardian REAL-Time monitor can be programmed to automatically receive fingerstick blood glucose calibration values transmitted from the LifeScan OneTouch® UltraLink™ Meter or the BD Paradigm Link® glucose meters (calibration values from these meters can also be manually entered). Calibrations from all from other glucose meters must be manually entered into the monitor.

The MiniLink transmitter is indicated for use as a component of select Medtronic Diabetes continuous glucose sensing systems and Medtronic Diabetes sensor-enabled pump systems.

Indications: CSS7100K (pediatric version)

The Guardian REAL-Time system (CSS7100K) is indicated for continuous or periodic monitoring of glucose levels in the fluid under the skin, in children and adolescents (ages 7 through 17) with diabetes mellitus, for the purpose of improving diabetes management. It alerts if a glucose level falls below, or rises above, pre-set values. Values are not intended to be used directly for making therapy adjustments, but rather to provide an indication of when a meter blood glucose measurement may be required. All therapy adjustments should be based on measurements obtained using a home glucose meter and not on Guardian REAL-Time CGM System values.

The Guardian REAL-Time CGM System (CSS7100K) provides real-time glucose values that allow users to track patterns in glucose concentrations and to possibly identify episodes of low and high blood glucose. It also stores the data so that it can be analyzed to track patterns. Glucose data can be further downloaded to PC software for analysis of historical glucose values. The Guardian REAL-Time monitor can be programmed to automatically receive

fingerstick blood glucose calibration values transmitted from the LifeScan OneTouch® UltraLink™ Meter or the BD Paradigm Link® glucose meters (calibration values from these meters can also be manually entered). Calibrations from all other glucose meters must be manually entered into the monitor.

The MiniLink transmitter is indicated for use as a component of select Medtronic Diabetes continuous glucose sensing systems and Medtronic Diabetes sensor-enabled pump systems.

Contraindications

Glucose monitoring therapy is not recommended for people who are unwilling or unable to perform a minimum of **two** meter blood glucose tests per day and to maintain contact with their healthcare professional. Successful monitor therapy requires sufficient vision or hearing to allow recognition of the monitor signals and alarms.

Safety Warnings

Sensor

The sensor can create special needs about your medical conditions or medications. Please discuss these conditions and medications with your healthcare professional before using the sensor.

Transmitter

Product contains small parts and may pose a choking hazard for young children.

The transmitter should be disconnected from the sensor while traveling on an aircraft or if it interferes with another transmitting device.

Meter

Do not use an RF-capable meter to send blood glucose measurements to your monitor while traveling on an aircraft or if it interferes with another transmitting device.

Precautions

Although the monitor has multiple safety alarms, it cannot tell you about all potential problems. You must test your blood glucose levels at least two times per day. If your blood glucose is out of range, then check your blood glucose using your blood glucose meter.

Avoid extreme temperatures

- Avoid exposing the system to temperatures above 108°F (42°C) or below 34°F (1°C).
- Do not steam, sterilize, or autoclave your monitor, transmitter, or sensor.

Sensor

Make sure the sensor is firmly attached before exercising or doing other physical activity.

Meter

Blood glucose measurements taken when your glucose levels change rapidly should not be used to calibrate your system. If your meter is set up to automatically transmit fingerstick blood glucose calibration values to your system, it is recommended that you turn off this feature before taking a fingerstick measurement when your blood glucose may be changing rapidly.

CAUTION: Any changes or modifications to the devices not expressly approved by Medtronic may damage the system and could void your warranty.

Adverse reactions

Operation of the sensor feature requires the insertion of a glucose sensor into the skin. Bleeding, swelling, bruising, or infection at the sensor insertion site are possible risks of sensor use. Contact your healthcare professional and your local representative in the event of any adverse reaction. The optional occlusive dressing should be removed if irritation from or reaction to the dressing develops. The sensor should be removed if redness, pain, tenderness or swelling develop at the insertion site.

Accuracy of the glucose sensor readings provided by the Guardian REAL-Time CGM System was not evaluated in patients with medical conditions, such as insulin resistance and cystic fibrosis. Also, the impact of medications was not specifically evaluated; therefore, device performance may vary in these situations.

Damage due to radiation

If you are going to have an X-ray, CT scan, MRI or be exposed to other types of radiation, REMOVE YOUR MONITOR, TRANSMITTER, and SENSOR and leave them outside the radiation area. The Guardian REAL-Time CGM System monitor and transmitter are designed and tested to tolerate common electromagnetic interference, including microwave ovens and airport security systems.

Compliance information

CAUTION: Changes or modifications not expressly approved by the party responsible for compliance could void the user's authority to operate the equipment.

The monitor, meter, and transmitter comply with the United States Federal Communications Commission (FCC) and international standards for Electromagnetic Compatibility. For the specific regulations and test results for your area, please contact your local representative.

The transmitter should be disconnected from the sensor while traveling on an aircraft, or if it interferes with another transmitting device. Do not use an RF-capable meter to send blood glucose measurements to your monitor while traveling on an aircraft. Manually enter your blood glucose reading into your monitor.

These devices comply with Part 15 of the FCC Rules. Operation is subject to the following two conditions:

- 1 These devices may not cause harmful interference
- 2 These devices must accept any interference received, including interference that may cause undesirable operation.

These standards are designed to provide reasonable protection against excessive radio frequency interference and prevent undesirable operation of the device from unwanted electromagnetic interference.

These devices generate, use, and can radiate radio frequency energy and may cause interference to radio or television communications, even if installed and used according to the instructions. If these devices do cause interference, try to correct the interference by one or more of the following measures:

- Reorient or relocate the monitor, meter and/or transmitter.
- Increase the separation between the monitor, meter, transmitter and the device that is receiving or emitting interference.

The transmitter and meter send information to the monitor using radio frequency. If other radio-frequency devices are in use, such as cell phones, cordless phones and wireless networks, they may prevent or interrupt communication between the monitor and the meter and/or the monitor and transmitter. This interference will not cause any incorrect data to be sent and will not cause any harm to your monitor, transmitter or meter. Moving away from or turning off these other devices will usually allow communication. See the *Troubleshooting and alarms/alerts* chapter to correct interference problems you may have. The transmitter does not interfere with any radio frequency signals transmitted from outside sources.

RF interference from other devices

Common consumer electronic devices that transmit in the same frequency band used by the MMT-7703 MiniLink transmitter may prevent the Guardian REAL-Time monitor from receiving the glucose information sent by the transmitter. Most cellular (mobile) phones and 900 MHz cordless phones, when transmitting or receiving, may cause significant interruption of transmitter-receiver communication. It is likely that other devices operating in similar frequency ranges will have a similar effect. This interference, however, will not cause any incorrect data to be sent and will not cause any harm to your transmitter

The Guardian REAL-Time monitor includes a programmable WEAK SIGNAL alert that notifies you when one or more expected transmissions were not received as expected by the receiving device. (The receiving device will also issue a LOST SENSOR alert if communication is interrupted for approximately 40 minutes.)

Communication problems can usually be resolved by:

- ensuring that there is less than 6 feet (1.8 meters) between the receiving device and transmitter, to include RF glucose meters;
- turning off or moving away from other RF transmitting devices; or
- reorienting or relocating the transmitter or receiving device, or both.

Communications problems may also occur due to cellular telephone interference. Testing conducted with several different cellular telephones indicates that, while a glucose value is being transmitted, using a cell phone within 12 inches (31 cm) of receiving devices, transmitters or RF glucose meters can interfere with reception of the transmitted values. In the event of such interference, normal communication can be re-established by:

- turning the cell phone off; or
- keeping the cell phone at least 12 inches away from the receiving device, transmitter or glucose meter when a glucose measurement is being transmitted.

If you have questions, please contact your local representative for assistance. Refer to the enclosed Medtronic Diabetes International Contact list for the HelpLine in your area.

The basics

In this chapter:

- The monitor
- Installing the battery
- The monitor buttons
- Alarms/alerts
- Operating modes
- Menus

CAUTION: Never use sharp objects to press the buttons on your monitor as this can damage the buttons or compromise the seal on the monitor. Some examples of sharp objects that can damage the buttons or the seal are fingernail files, pens, pencils, paper clips, knives, scissors, and keys.

The monitor

The screen shows information such as your real-time sensor glucose measurements, alarm/alert messages, current settings, and status readings. The buttons are used to select menu items, enter values, and perform actions.



Installing the battery

CAUTION: Do not use a rechargeable or carbon zinc battery in your monitor. For best results, a new AAA Energizer alkaline battery is recommended.

Medtronic has designed an important safety measure into the monitor: It measures battery strength during installation. If you install a new, full-strength battery, the monitor accepts it. However, if you install a battery that does not have full power, or take too long to switch batteries, the monitor will show one of the following alarms:


- **WEAK BATTERY** – The battery is weak but usable. Follow the instructions on the screen to clear the alarm, then either install the weak battery, or install another new battery. Your monitor will work with the weak battery, but the battery life will be reduced.
- **FAILED BATT TEST** – The battery is not strong enough to power your monitor; your monitor will not work with that battery. Clear the alarm and install a new battery.
- **BATT OUT LIMIT** – May occur if the battery is removed from the monitor, and more than 10 minutes passes before the new battery is inserted. If this alarm occurs, continue installing the new battery.

Battery life varies widely depending on the way you use your monitor. Some of the use factors that affect battery life include:

- the alert type you choose
- the number of alerts you receive
- how much you use the backlight
- how long you display the graphs


To install the battery:

- 1 If your new battery has been stored in a refrigerator or other cold environment, you must let it warm up to room temperature before installing it in your monitor.
- 2 Before you remove the old battery, make sure that you set the monitor on the HOME (idle) screen and clear any alarms (ESC, ACT).
- 3 Use the edge of a coin to remove the battery cap. Turn the cap counter-clockwise.

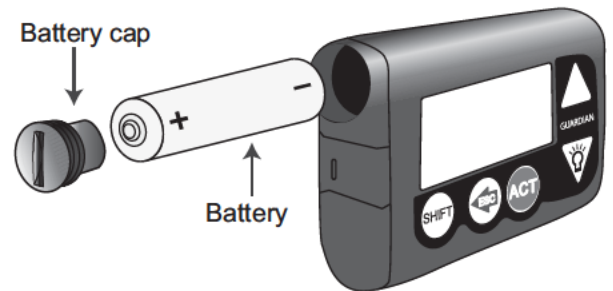


- 4 Remove the old battery and throw it away according to the laws of your state or country. Insert the new battery into the monitor negative end [(-) symbol] first. Check the label on the back of the monitor to make sure that the battery is inserted correctly.
- 5 Place the battery cap in the monitor and tighten until the slot is aligned with the monitor as shown at right.

CAUTION: Do not over-tighten the battery cap. You should not turn the cap more than four half-turns. If you over-tighten the cap you may not be able to remove it, and you can damage your monitor.

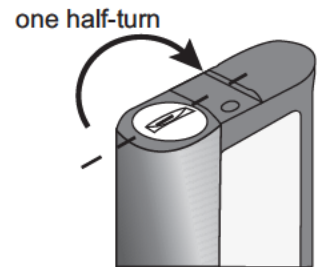

When the monitor turns on, several screens appear, the monitor shows the version number, then the HOME screen appears. If you are replacing a battery for a monitor that has very recently been in use, the only screen you will see before the HOME screen is the one showing the version number.

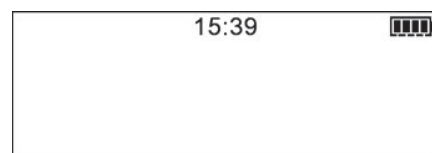

- 6 If the HOME screen appears, then go to the next step. If the HOME screen does not appear:
 - a. Make sure that the battery is inserted correctly. If not, remove the battery and install it properly as described in the previous steps.
 - b. If the monitor still does not turn on, or if you receive a FAILED BATT TEST alarm, then remove the battery and replace it with a new one.
 - c. If the monitor still does not turn on, then call the HelpLine.
- 7 Make sure that the time showing on the HOME screen is correct. If more than 10 minutes have passed since you removed the battery, you will be prompted to check the time and date.
- 8 Press **ESC** until the STATUS screen appears. Make sure that no alarms or alerts are active. If an alarm or alert is active, then follow the instructions on the screen to clear the alarm or alert.

The monitor buttons

This section describes how to use the buttons on your monitor.

Button	Description
	• Scroll up through menus when selecting menu items. • Increase the value of a setting when entering or changing settings. • Select different graphs when viewing the sensor glucose graphs. • Move the cursor to the right in the graph when viewing the sensor glucose graphs.
	• Scroll down through menus when selecting menu items. • Decrease the value of a setting when entering or changing settings. • Move the cursor to the left in the graph when viewing the sensor glucose graphs. • Turn the backlight on when the HOME screen is showing. • Turn the backlight on when viewing screens other than the HOME screen (hold down SHIFT and press DOWN).
	• Accept the selected menu item. • Choose a setting when entering or changing settings. • Open the MAIN MENU from the HOME screen. • Clear alarms/alerts (press ESC, then press ACT).
	• Return to previous screen or menu. • Open the following screens and graphs from the HOME screen: • The sensor glucose graphs. • The STATUS screen and SENSOR STATUS screen. • The SENSOR DEMO screens (if the feature is turned On). • Move the cursor to the far right side of the graph when viewing the sensor glucose graphs. • Clear alarms/alerts (press ESC, then press ACT).
	• Turn the backlight on when viewing screens other than the HOME screen (hold down SHIFT and press DOWN). • Access the User Settings (hold down SHIFT and press ACT).

An important message about alarms and alerts

Your Guardian REAL-Time system tells you about events and conditions that affect your diabetes management in two ways:

- With **alarms** about events and/or conditions that affect the ability of your system to work properly. An alarm means you must do something immediately to make sure your system keeps working properly. When an alarm occurs, a closed circle icon appears at the top left of your monitor screen.  closed circle
- With **alerts** about your glucose levels or about conditions that affect the way your system measures your glucose levels. An alert means your glucose levels are changing or you must do something within a certain amount of time to make sure your system's glucose monitoring stays accurate. When an alert occurs, an open circle icon appears at the top left of your monitor screen.  open circle


Certain alerts are called *alarms* on your screen. When this happens, the most important information is whether the circle icon on your screen is open (alert) or closed (alarm). See the *Troubleshooting and alarms/alerts* chapter for details on alarms and alerts.

The monitor screen

The screen shows five lines of text at one time. The first is the operating mode. The second is the current open menu or function. The last three lines show either information or text that you can select for the current function.

The HOME screen

Start from the HOME screen to get to the menus and screens you use to set up the system. Use the HOME screen to view the time and the sensor icon, which tells you if the transmitter and monitor are communicating. If you do not press any buttons for approximately 30 seconds, the monitor automatically returns to the HOME screen.

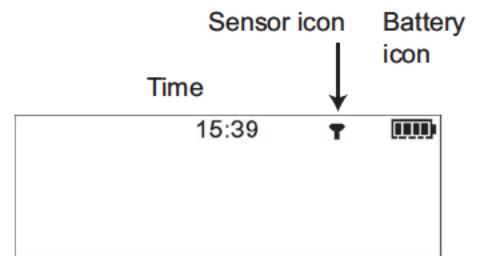

When the monitor is on, the time and the battery icons always show across the top of the HOME screen. The sensor icon will not show until you start the sensor and establish communication between the monitor and transmitter. If the time and battery icons do not show, the monitor is not working.

Icons

The icons that show on the monitor screen are described below.

Sensor icon

A dark sensor icon means that the sensor is on and the monitor is receiving sensor data. A light-colored sensor icon with a dark circle around it means that the sensor is on, but the monitor is not receiving sensor data. For information about the sensor and transmitter, see the *Setting up your sensor, transmitter and meter* chapter.




Sensor is on and the monitor is receiving sensor data



Sensor is on but the monitor is not receiving sensor data

Battery icon

This icon shows when a battery is installed in the monitor. There are four segments in the icon. Each segment represents approximately 25 percent of the usable battery life you have left until you reach Low Battery point. If you only have one segment left, make sure you have a new battery available.

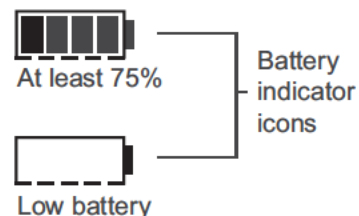


Time of day icon

The current time of day shows across the top of the monitor screen in the format you select — 12-hour or 24-hour. See the *Entering your settings* chapter for details on setting the time.


Alert icon

An open circle icon appears when an alert occurs. See the *Troubleshooting and alarms/alerts* chapter for information on the different alerts.

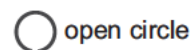

Alarm icon

A closed circle icon appears when an alarm occurs. See the *Troubleshooting and alarms/alerts* chapter for information on the different alarms.

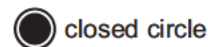

Scroll bar

A scroll bar appears on the right side of the screen when there is more information on the screen than the monitor can show at the same time. Press  to scroll to the other information.

Screen backlight

Press  to turn the screen backlight on or off in the HOME screen. If you are in another screen, hold **SHIFT**, then press  to turn on the backlight. The light will stay on as long as you use the current screen.

To conserve your battery, the backlight turns off automatically when the monitor vibrates and turns on again when the vibration ends. The backlight will not turn on in a Low Battery condition.

Alarm/alert types


You can set up your monitor to either beep or vibrate when an alarm or an alert occurs. See the next chapter for details on setting the alert type.

Operating modes

The screen tells you when a special feature is active or if there is a condition that needs your attention. The monitor status will determine the operating mode. The screens for the different modes are shown below.

Normal mode


No alarms or alerts are active.

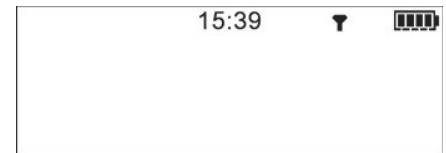

Special mode (Alerts)


Tells you that a special feature is active or that a certain alert condition exists. In Special mode, an open circle appears at the top of the screen and the monitor beeps or vibrates periodically to remind you of the condition. All sensor alerts change from beep or vibration to siren if not cleared within 10 minutes. Special mode does not restrict any monitor functions. Special mode features and conditions are:

- Block feature is turned on (special feature)
- Low Battery alert
- Sensor alerts

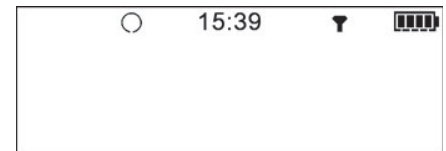

Attention mode (Alarms)

Tells you that an alarm is active or that there is an alarm condition that needs immediate attention. In the Attention mode, a solid circle shows at the top of the screen and the monitor beeps or vibrates periodically until the condition is cleared. The beep or vibration frequency changes, depending on the condition that put the monitor in Attention mode.

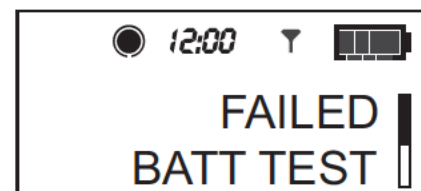

The screen shows the name of the alarm and information about how to clear it (press  to see the information on clearing the alarm). In this example, the alarm is the FAILED BATT TEST alarm.

Menus

There are three menus. Each menu contains different menu items that you select to set up the system and to view information.

MAIN MENU

This is the first menu in the system. It contains the following sub-menus: SENSOR, CAPTURE EVENT, and UTILITIES. To open the MAIN MENU, press ACT from the HOME screen.

SENSOR MENU

Contains the options that you use to set up communication between the sensor and the monitor, and to get to the monitor's sensor features. See the next chapter for details on how to enter your glucose monitoring settings. See the *Setting up your sensor, transmitter and meter* chapter for details on how to set up your sensor.

CAPTURE EVENT

Contains the options that you use to enter information about different events into the system. Examples of events include: the amount of insulin used for injections; and carbohydrates consumed during meals or snacks. See the *Using your system* chapter for details on how to enter event information.

UTILITIES MENU

Contains the options that you use for safety features, convenience, and to view system alarm/alert information. Examples include: locking the keypad buttons; viewing system alarms/alerts; turning the meter option on and off; and running the Selftest diagnostic test. See the *Utilities* chapter for details on using this menu.

The STATUS screen

The STATUS screen shows any recent activity, the condition of the monitor battery, and which special features are active, if any.

- To open the STATUS screen from the HOME screen, press **ESC** until the STATUS screen appears.
- To view additional information, press .
- To return to the HOME screen, press **ESC** until the HOME screen appears.

For details on the information shown in the STATUS screen, see the *Using your system* chapter.

Do NOT try to open the STATUS screen when you enter settings into the monitor. If you press the **ESC** button when you are entering settings, you may cancel the settings.

Entering your settings

In this chapter:

- Basic steps
- Setting the time and date
- Selecting the language
- Setting your Alarm/Alert type
- Entering your glucose monitoring settings
- Reviewing your settings

Basic steps

Perform these basic steps in the order shown.

- 1 Make sure that you have installed the battery in your monitor.
- 2 Make sure that you have the information you need to begin.
- 3 Set the time and date.
- 4 Set the alert type.
- 5 Enter your glucose monitoring settings.
- 6 Review your settings.
- 7 Save your settings.
- 8 Modify your settings as needed.

Setting the time and date

You must set the correct time and date in your monitor so that the system accurately reflects the time and date of your reported glucose values. You must reset the time and date if you receive a CHECK SETTINGS alarm, clear your settings, or remove the battery for more than ten minutes. You can select a 12-hour or 24-hour clock. If you do not read English, you must set the time and date using the English menus before you can select the language. For your convenience, the English screens you use to set the time and date are labeled in your language.

- 1 From the HOME screen, open the TIME/DATE SETUP screen.

MAIN MENU > Utilities > Time/Date

The TIME/DATE SETUP screen appears.




- 2 Select **12-Hour Setup** or **24-Hour Setup**. Press **ACT**. The TIME/DATE SET screen shows your monitor's current time/date setting.
- 3 Press **ACT**. The SET HOUR screen appears.

4 Use  /  to set the correct time and date as follows:

- a. Set the correct hour, then press **ACT**.
- b. Set the correct minutes, then press **ACT**.
- c. Set the correct year, then press **ACT**.
- d. Set the correct month, then press **ACT**.
- e. Set the correct day, then press **ACT**.



- 5 If you are setting up your monitor for the first time, then the TIME/DATE CORRECT screen appears.
- a. Make sure that the time and date shown on the screen is correct.
 - b. If the time and date is **correct**, select **Yes**, then press **ACT**. If it is **not correct**, press **ACT**. The TIME/DATE SETUP screen reappears. Repeat the procedure and enter the correct information.
- 6 If you are re-setting the time only, then the TIME SET AT screen appears.
- a. Make sure that the time and date shown on the screen is correct.
 - b. If the time and date is correct, press **ACT**. If it is not correct, press **ESC** until the HOME screen appears, then repeat the procedure to enter the correct time and date.

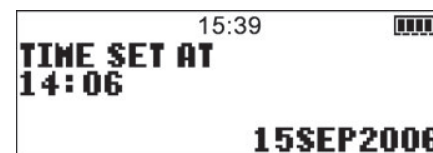


Selecting the language

If you do not read English, the next step is to change the Language setting so that all text information, alarms, alerts, sensor updates, meter blood glucose measurements, and sensor glucose measurements show on your screen in your language (some languages are not available on all monitors). You must make this change using the English menus until the language setting is actually changed.

To change the language setting:

- 1 Open the LANGUAGE MENU.
MAIN MENU > Utilities > Language
- 2 Press **ACT**. The screen shows English selected as the default language.
- 3 Scroll up or down to select the language you want to use.
- 4 Press **ACT**. The language setting is changed and the screen returns to the bottom of the UTILITIES menu in the selected language.
- 5 To return to the HOME screen, press **ESC** until it appears.

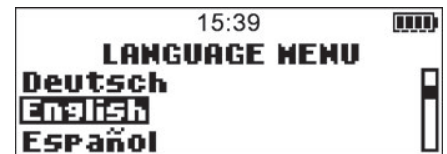

Setting your Alarm/Alert Type

You can set your monitor three different ways to let you know when there is an alarm, alert, or special conditions or programming (see the *Troubleshooting and alarms/alerts* chapter for details). The default setting is Beep Medium, but you can select Beep Long, Beep Short, or Vibrate.

The Vibrate setting uses more battery power than beep settings and may shorten battery life. If your alarm/alert type is set to Vibrate and you get a LOW BATTERY alert, then your monitor will automatically use the default beep to save battery power.

To set the alert type:

- 1 Open the ALERT TYPE screen.

MAIN MENU > Utilities > Alarm > Alert Type

The ALERT TYPE screen shows with Beep Medium selected (the default alert type).

- 2 Select your alert type (Beep Short, Beep Medium, Beep Long, Vibrate). You must scroll down to select Vibrate.
- 3 Press **ACT**. The alert type is set. Alarms/alerts will beep or vibrate as selected.
- 4 Press **ESC** until the HOME screen appears.

Entering your glucose monitoring settings

Enter your glucose monitor settings in the exact order they appear below, as some of the settings are dependent upon other settings being made first. When you complete each setting, your monitor will automatically display the screen for the next setting in the sequence. When a screen item flashes, you may change the value of the flashing item by pressing  or .

Turning on the sensor

The sensor must be turned on, started and initialized to report glucose measurements.

- 1 Open the EDIT SETTINGS screen.

MAIN MENU > Sensor > Sensor Setup > Edit Settings

The EDIT SETTINGS screen appears with **Sensor: Off** selected.

- 2 Press **ACT**. The SENSOR ON/OFF screen appears with **Off** selected.
- 3 Select **On**, then press **ACT**. The EDIT SETTINGS screen appears showing the sensor is now turned on.

Selecting the BG units

You must select the blood glucose units (BG units) for the system (either mg/dL or mmol/L). All BG measurements will show in the BG unit type you select.

- 1 Make sure the EDIT SETTINGS screen is open.
MAIN MENU > Sensor > Sensor Setup > Edit Settings
- 2 Select **BG Units**, then press **ACT**. The BLOOD GLUCOSE UNITS screen shows **mg/dL** or **mmol/L**.
- 3 Press  or  to select the units you want.
- 4 Press **ACT**. The EDIT SETTINGS screen shows the BG Units type you selected.

Turning on the Glucose Alerts

The Glucose alerts must be turned on so the system can send you an alert when the sensor glucose measurements reach or exceed your Glucose Limits.

- 1 Make sure the EDIT SETTINGS screen is open.
MAIN MENU > Sensor > Sensor Setup > Edit Settings
- 2 Select **Glucose Alarms**, then press **ACT**. The GLUCOSE ALARMS screen appears with **Off** selected.
- 3 Select **On**, then press **ACT**. The EDIT SETTINGS screen appears showing that the glucose alarms are now turned **On**.

Setting the Glucose Limits


After you turn the Glucose alerts **On**, then you need to set the high and low Glucose Limits recommended by your healthcare professional. Your recommended glucose limits may vary throughout the day, so your monitor allows you to set up to eight pairs.

A Glucose Limit pair includes one High Glucose Limit and one Low Glucose Limit, as shown in the following table:

Glucose Limit pair	Alerts
Low: 80 mg/dL (4.4 mmol/L) for CSS7100 90 mg/dL (5.0 mmol/L) for CSS7100K High: 200 mg/dL (11.1 mmol/L) for CSS7100 280 mg/dL (15.6 mmol/L) for CSS7100K	• A Low Glucose Limit alert occurs when the sensor glucose measurement reaches or goes below the Low Glucose Limit in this table. • A High Glucose Limit alert occurs when the sensor glucose measurement reaches or goes above the High Glucose Limit in this table.

Setting the Glucose Limits start time

If your healthcare professional recommends that you set more than one pair of Glucose Limits, you must enter a start time for each. For example, you may use one pair of Glucose Limits during the day and another pair at night. The example below shows two pairs of Glucose Limits. The first pair begins at midnight, and the second pair begins at 9:00.

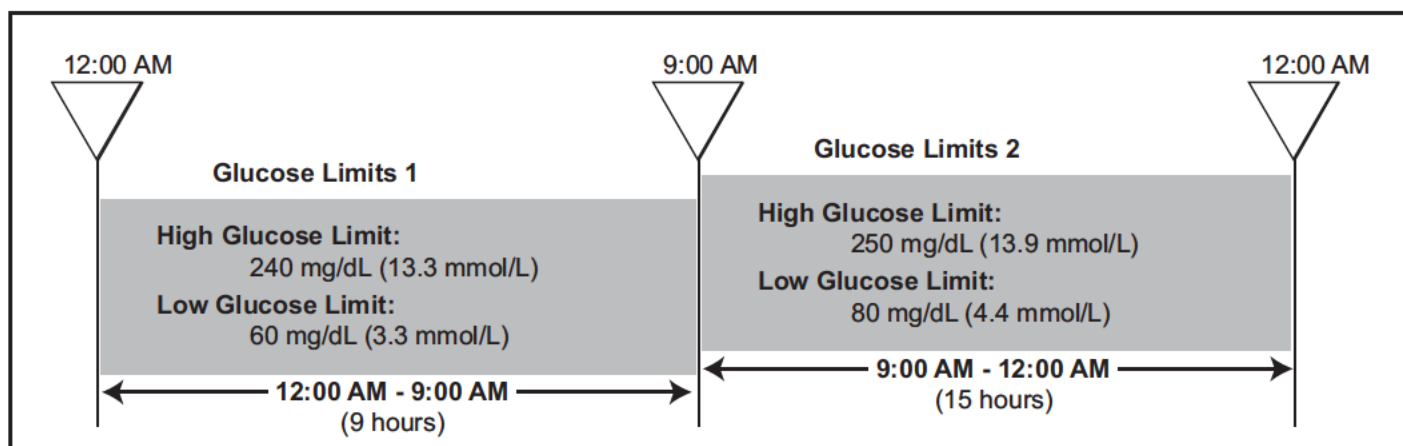

To set up your Glucose Limits:

- 1 Make sure the EDIT SETTINGS screen is open.

MAIN MENU > Sensor > Sensor Setup > Edit Settings

- 2 Select Glucose Limits, then press ACT. The SET GLUCOSE LIMITS 1 screen appears. The Low Glucose Limit flashes:
 - 60 mg/dL (3.3 mmol/L) for CSS7100
 - 90 mg/dL (5.0 mmol/L) for CSS7100K

The start time for your first pair of Glucose Limits is midnight (00:00) and cannot be changed.

- 3 Select your Low Glucose Limit. The available range is 40 to 390 mg/dL (2.2 to 21.6 mmol/L).
- 4 Press ACT. The High Glucose Limit flashes:
 - 200 mg/dL (11.1 mmol/L) for CSS7100
 - 280 mg/dL (15.6 mmol/L) for CSS7100K
- 5 Select your High Glucose Limit. The available range is 50 to 400 mg/dL (2.8 to 22.2 mmol/L). Your High Glucose Limit must be at least 10 mg/dL (0.6 mmol/L) above your Low Glucose Limit.

30 Entering your settings

- 6 Press **ACT**. The SET START TIME 2 screen appears. The dashes under the screen name flash. The first pair of Glucose Limits is now set.
- 7 To set a second pair of Glucose Limits, select the time of day you want this pair to become active. For example, if you select 9:00, then this pair of Glucose Limits will become active at 9:00.
- 8 Press **ACT**. The SET GLUCOSE LIMITS 2 screen appears.
- 9 The default Low Glucose Limit flashes (**OFF**).
- 10 Select your Low Glucose Limit. The available range is:
 - 40 to 390 mg/dL (2.2 to 21.6 mmol/L) for CSS7100
 - 90 to 390 mg/dL (5.0 to 21.6 mmol/L) for CSS7100K
- 11 Press **ACT**. The default High Glucose Limit flashes (**OFF**).
- 12 Select your High Glucose Limit, then press **ACT**. The value must be between:
 - 50 and 400 mg/dL (2.8 and 22.2 mmol/L) for CSS7100
 - 100 and 400 mg/dL (5.56 and 22.2 mmol/L) for CSS7100K
- 13 Press **ACT**. The SET START TIME 3 screen appears.
- 14 The second pair of Glucose Limits (Glucose Limits 2) is now set. If you do not need any more Glucose Limits pairs, press **ESC**. If you want to set more Glucose Limits, then repeat step 7 through step 13 above for each pair, up to a maximum of eight.
- 15 If you wish to switch off either the Low or High Glucose Limits independently of each other, reduce the Low Glucose Limit below 40 mg/dL (2.2 mmol/L) for step 3 to select **OFF**, or reduce the High Glucose Limit below 50 mg/dL (2.7 mmol/L) for step 5 to select **OFF**.

Setting the High Snooze

After you receive and clear a High Glucose, Rise Rate of Change, or High Predictive alert, the alert will repeat until the condition that caused the alert is resolved. The High Snooze feature allows you to set how frequently you want the alert to repeat after you clear it the first time. For example, if your glucose levels usually take one to two hours to fall after a correction bolus or injection, you may want to set the High Snooze to one or two hours to avoid unnecessary alerts.

Example

You set your High Snooze to 15 minutes. You receive an alert at 13:00 and immediately clear it; however, the condition that caused the alert still exists, so you will receive a second alert at 13:15, 15 minutes after you cleared the first alert. If you immediately clear the alert again, then the alert will repeat at 13:30 and continue to repeat every 15 minutes until you resolve the condition that caused the alert. (This example also applies to the Low Snooze.)

To set up your High Snooze:

- 1 Make sure the EDIT SETTINGS screen is open.
MAIN MENU > Sensor > Sensor Setup > Edit Settings
- 2 Select **High Snooze**, then press **ACT**. The SET HIGH SNOOZE screen flashes the default High Snooze time of 1 hour (1:00).
- 3 Select your High Snooze time. The time must be between five minutes (0:05) and three hours (3:00).
- 4 Press **ACT**. The EDIT SETTINGS screen shows the High Snooze time set.

Setting the Low Snooze

After you receive and clear a Low Glucose, Fall Rate of Change, or Low Predictive alert, the alert will repeat until the condition that caused the alert is resolved. The Low Snooze feature allows you to set how frequently you want the alert to repeat after you clear it the first time.

- 1 Make sure the EDIT SETTINGS screen is open.
MAIN MENU > Sensor > Sensor Setup > Edit Settings
- 2 Select **Low Snooze**, then press **ACT**. The SET LOW SNOOZE screen flashes the default Low Snooze time of 20 minutes (0:20).
- 3 Select your Low Snooze time. The time must be between 5 minutes (0:05) and 1 hour (1:00).
- 4 Press **ACT**. The EDIT SETTINGS screen shows the Low Snooze is set.

Setting the Predictive glucose alert

The Predictive alerts calculate when you will reach your Low or High Glucose Limits, then send you an alert before you reach those limits. A Predictive alert tells you that **if your sensor glucose measurements keep falling or rising at the current rate**, you will reach your Glucose Limit in the number of minutes you selected.

Understanding the Time Sensitivity Settings

Time Sensitivity means the length of time, in minutes, you want to be notified before you reach your High or Low Glucose Limit. If you set the Low time sensitivity setting at 25 minutes and the High time sensitivity setting at 20 minutes, then the Predictive alerts will be sent:

- 25 minutes **before** the sensor glucose reaches your Low Glucose Limit
- 20 minutes **before** the sensor glucose reaches your High Glucose Limit

To set up the Predictive alerts:

- 1 Make sure the EDIT SETTINGS screen is open.
MAIN MENU > Sensor > Sensor Setup > Edit Settings
- 2 Select **Predictive Alm**, then press **ACT**. The PREDICTIVE ALARMS LOW/HIGH screen shows that the alerts are turned **Off**.
- 3 Select **On**, then press **ACT**. The SET TIME SENSITIVITY LOW/HIGH screen appears and the LOW time setting flashes the default time: 15 minutes.
- 4 Select the Time Sensitivity (how many minutes before you reach your Low Glucose Limit).
 - Range: **5 to 30 minutes (0:05 to 0:30)** 5 to 30 minutes (0:05 to 0:30)
 - Increments: 5 minutes
- 5 Press **ACT**. The HIGH time setting flashes the default time: 15 minutes.
- 6 Select the Time Sensitivity (how many minutes before you reach your High Glucose Limit).
 - Range: **5 to 30 minutes (0:05 to 0:30)** 5 to 30 minutes (0:05 to 0:30)
 - Increments: 5 minutes
- 7 Press **ACT**. The EDIT SETTINGS screen appears.

Setting the Rate Alerts

The Rate alerts tell you when your sensor glucose (SG) changes at, or faster than, the per-minute rate pre-selected by you. There are two alerts:

- **FALL RATE** for SG **decreases** at or faster than your pre-selected rate
- **RISE RATE** for SG **increases** at or faster than your pre-selected rate

You can set one or both Rate alerts. The rate alert default setting is **OFF**. When you press  or , the number 4.0 mg/dL (0.220 mmol/L) shows, then adjust up or down with each additional press.

The explanation below shows how the Rate alerts work.

If you set the number	Then
Lower than 4.0 mg/dL (0.220 mmol/L)	• The system is more sensitive to SG changes • Alerts are more frequent than if you use a higher number
Higher than 4.0 mg/dL (0.220 mmol/L)	• The system is less sensitive to SG changes • Alerts are less frequent than if you use a lower number

To set up the Rate alerts:

- 1 Make sure the EDIT SETTINGS screen is open.
MAIN MENU > Sensor > Sensor Setup > Edit Settings
- 2 Select **Rate Alarms**, then press **ACT**. The SET FALL RATE LIMIT screen flashes the rate alert default setting: **OFF**.
- 3 Select the Fall Rate Limit.
 - Range: mmol/L/min from **0.061** to **0.277**
 mg/dL/min from **1.1** to **5.0**
 - Increments: mmol/L/min **0.005**
 mg/dL/min **0.1**
- 4 Press **ACT**. The SET RISE RATE LIMIT screen flashes the rate alert default setting: **OFF**.
- 5 Select the Rise Rate Limit. The range and increments are the same as the Fall Rate Limit (see step 3 above).
- 6 Press **ACT**. The Rate alerts will now use your settings.

Setting the AUC Limits

The AUC (Area Under the Curve) Limits feature measures how much and how long your sensor glucose measurements are above the pre-defined AUC High Limit or below the pre-defined AUC Low Limit.

Check with your healthcare professional to determine the AUC settings that are best for you. The AUC Limits are used to analyze the sensor glucose measurement data that is stored in your system. The AUC Limits are not high and low glucose alerts. AUC Limits are different than Glucose Limits.

The default AUC Limits are:

- Low AUC Limit: 70 mg/dL (3.9 mmol/L)
- High AUC Limit: 180 mg/dL (10 mmol/L)

To set the AUC Limits:

- 1 Make sure the EDIT SETTINGS screen is open.
MAIN MENU > Sensor > Sensor Setup > Edit Settings
- 2 Select **AUC Limits**, then press **ACT**. The SET AUC LIMITS screen appears. The default LOW AUC Limit flashes 70 mg/dL (3.9 mmol/L).
- 3 Select your Low AUC Limit. The number must be between 40 and 400 mg/dL (2.2 and 22.2 mmol/L).
- 4 Press **ACT**. The default High AUC Limit flashes 180 mg/dL (10.0 mmol/L).
- 5 Select your High AUC Limit. The number must be between 40 and 400 mg/dL (2.2 and 22.2 mmol/L). It can be equal to, but not lower than, the AUC Low Limit.
- 6 Press **ACT**. The EDIT SETTINGS screen shows the AUC Limits you selected.

Setting the Alarm Snooze

After you receive and clear a Meter BG Now alert, the monitor will repeat the alert until you enter a new meter blood glucose measurement. The Alarm Snooze feature allows you to set how frequently you want the alert to repeat after you clear it.

- 1 Make sure the EDIT SETTINGS screen is open.
MAIN MENU > Sensor > Sensor Setup > Edit Settings
- 2 Select **Alarm Snooze**, then press **ACT**. The SET ALARM SNOOZE screen flashes the default Alarm Snooze time of 30 minutes (0:30).
- 3 Set your Alarm Snooze time between 5 minutes (0:05) and 1 hour (1:00).
- 4 Press **ACT**. The EDIT SETTINGS screen shows the Alarm Snooze time you selected.

Setting the Cal Reminder

After you calibrate your system for the first time, you must calibrate it again within the next six hours. The Cal Reminder feature allows you to set a reminder to calibrate your system. For example, if you set your reminder to four hours, then you will receive a Cal Reminder (Meter BG By) four hours before the next meter BG entry is due (eight hours after your last successful sensor calibration). If you miss a calibration, the system will send you a METER BG NOW alert.

To set up the Cal Reminder:

- 1 Make sure the EDIT SETTINGS screen is open.
MAIN MENU > Sensor > Sensor Setup > Edit Settings
- 2 Select **Cal Reminder**, then press **ACT**. The CAL REMINDER screen shows **On**.
- 3 Press **ACT**. The SET CAL REMINDER screen flashes the Cal Reminder default time of 1 hour (1:00).
- 4 Select your Cal Reminder time between 5 minutes (0:05) and 6 hours (6:00).
- 5 Press **ACT**. The EDIT SETTINGS screen shows the Cal Reminder time you selected. You are now ready to enter the ID number of your transmitter.


Entering the transmitter identification number

The Transmitter ID starts with 2 and is found on the flat side of your transmitter. You must enter the transmitter identification number so that the transmitter and monitor can communicate with each other.

Write the Transmitter ID of the transmitter you are using here: _____.

- 1 Make sure the EDIT SETTINGS screen is open.
MAIN MENU > Sensor > Sensor Setup > Edit Settings
- 2 Select Transmitter ID, then press **ACT**. The SET TRANSMITTER ID screen appears. The screen shows seven dashes for the Transmitter ID, the first of which (on the left) is flashing.
- 3 Select the first number of the Transmitter ID (the number on the screen must match the first number of the Transmitter ID). The transmitter ID is on the transmitter label, as shown.
- 4 Press **ACT**. The second dash flashes.

SN GT2XXXXXXM

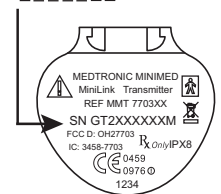

- 5 Select the next number of the Transmitter ID, then press **ACT**.
- 6 Repeat step 5 until all seven digits of the Transmitter ID have been selected.
A message shows that the Transmitter ID has been changed. After approximately 30 seconds, the HOME screen appears.

Setting up the Missed Data option

The Missed Data option allows you to determine how quickly you will receive a Weak Signal alert when there is no communication between the transmitter and the monitor. For example, if you set the Missed Data alert to 30 minutes, you will receive a Weak Signal alert when your transmitter and monitor have been out of communication for 30 minutes. The transmitter has 40 minutes of memory.

To set up the Missed Data option:

- 1 Make sure the EDIT SETTINGS screen is open.
MAIN MENU > Sensor > Sensor Setup > Edit Settings
- 2 Select **Missed Data**, then press **ACT**. The SET MISSED DATA screen flashes the default time of 30 minutes (0:30).
- 3 Select your Missed Data time between 5 minutes (0:05) and 40 minutes (0:40).
- 4 Press **ACT**. The EDIT SETTINGS screen shows the length of time you entered for the Missed Data option.

Setting the sensor glucose graph timeout

The sensor glucose graph timeout is how long that graphs show on the monitor screen until the HOME screen automatically appears. The default timeout is two minutes, and you can set it for 2, 4, or 6 minutes.

Using a continuous graph display

You can also select NONE to set the sensor graph timeout to display continuously. If you select NONE, the sensor graphs will show until another screen or menu is selected, or an alarm/alert occurs. Selecting continuous graph display uses more battery power and decreases battery life. To return to the HOME screen from the sensor glucose graphs, press **ESC** until the HOME screen appears.

To set the sensor glucose graph timeout:

- 1 Make sure the EDIT SETTINGS screen is open.
MAIN MENU > Sensor > Sensor Setup > Edit Settings
- 2 Select **Graph Timeout**, press **ACT**. The SET GRAPH TIMEOUT screen flashes the default time (two minutes).
- 3 Select your Graph Timeout. The time you select must be: two minutes; four minutes; six minutes; or NONE.
- 4 Press **ACT**. The EDIT SETTINGS screen shows the Graph Timeout you selected.

Reviewing your settings

Make sure that all the settings you have made are correct before you use the system. Use the following procedure to review your settings:

- 1 Open the REVIEW SETTINGS screen.
MAIN MENU > Sensor > Sensor Setup > Review Settings
- 2 Press  to scroll through all settings to make sure that they are correct.
- 3 To change any settings, return to the EDIT SETTINGS menu.
MAIN MENU > Sensor > Sensor Setup > Edit Settings
- 4 Save your settings after you have set your preferences so that you can restore the settings you have saved if you receive an alarm/alert or error that resets your settings.

Using the Sensor Demo

The Sensor Demo shows you demonstrations of the sensor graphs and the Sensor Status screen. To view the Sensor Demo screens, you must first turn on the Sensor Demo feature.

Turning on the Sensor Demo feature

- 1 Make sure the EDIT SETTINGS screen is open.
MAIN MENU > Sensor > Sensor Setup > Edit Settings
The EDIT SETTINGS screen shows **Sensor Demo: Off**.
- 2 Press **ACT**. The SENSOR DEMO ON/OFF screen shows with **Off** selected.
- 3 Select **On**, then press **ACT**. The EDIT SETTINGS screen shows that SENSOR DEMO is now turned **on**.

Viewing the Sensor Demo screens

- 1 From the HOME screen, press **ESC**. The screen shows SENSOR DEMO, then the first Sensor Demo graph appears.
- 2 You can move the cursor (flashing vertical line) through the graph to see examples of how the actual real-time graphs will look at specific times. Press  to move the cursor left, press  to move the cursor right.
- 3 You must turn the Sensor Demo **Off** in order to show your actual, real-time sensor glucose graphs. To turn the Sensor Demo feature off, follow the above procedure and select **Off**.

Setting up your transmitter, sensor and meter

In this chapter:

- Basic steps
- Setting up the transmitter
- Inserting the sensor
- Connecting the transmitter to the sensor
- Starting the sensor
- Entering the first meter BG

Basic steps

After you have set up your monitor, there are a few basic steps that you must complete to set up your sensor and transmitter for use. Perform these steps in the order below.

- 1 Install the charger battery.
- 2 Charge the transmitter battery.
- 3 Insert the sensor into your body.
- 4 Connect the transmitter to the sensor.
- 5 Start the sensor.
- 6 Enter your first meter BG.

Setting up the transmitter


The MiniLink transmitter is a device that takes electronic signals generated by the glucose sensor and sends them by radio frequency to the monitor. It has a blue tester and a charger.

MiniLink charger

The MiniLink transmitter contains a non-replaceable, rechargeable battery that you can recharge as needed with the charger. The charger has a green light that shows the charging status and a red light that communicates any problems during charging. If you see a red light, see the *Troubleshooting and alarms/alerts* chapter. The charger needs one AAA battery to operate (Energizer is recommended). A new AAA battery contains enough power to recharge the MiniLink transmitter more than 40 times.

Installing a new charger battery

- 1 Remove the charger battery cap by turning it counter-clockwise 1/4 turn using a coin in the groove of the cap.

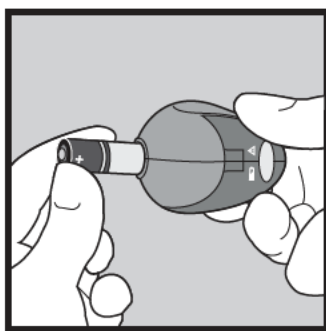


- 2 Insert a new AAA battery with the flat (-) end first. Make sure that you align the small bumps on the battery cap with the small notches in the charger's battery opening. Push the cap in all the way using a coin. Turn the cap clockwise 1/4 turn.
- 3 If the battery is installed incorrectly or is low, the charger will not work. Repeat the steps above using a new battery.


Charging the transmitter

Before using the transmitter for the first time, you must fully charge the transmitter battery, which may take up to eight hours. It is recommended to recharge the transmitter after each sensor use. If you choose to recharge the transmitter after each sensor use, the charging time will be less than 20 minutes. A fully charged transmitter battery will work more than 14 days without recharging. After 14 days of use, the transmitter will fully recharge in less than two hours.


- 1 If a green light on the transmitter is lit or flashing, do not connect it to the charger. The transmitter will not charge with its green light on. Wait for the green light to turn off, then connect the transmitter to the charger.
- 2 Connect the transmitter to the charger by lining it up, flat side down, with the charger. Push the two components together fully. Always allow at least one minute before disconnecting the transmitter from the charger or the transmitter may not work properly. If you disconnect the transmitter before one minute, reconnect it to the charger for at least one minute.

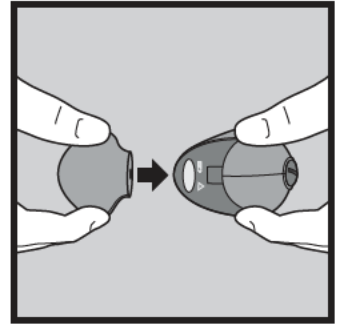

- 3 Within 10 seconds after the transmitter is connected, a green light on the charger will flash for one to two seconds as the charger powers on. For the rest of the charging time, the charger's green light will flash in a continuing pattern of four flashes, pause, four flashes, pause.

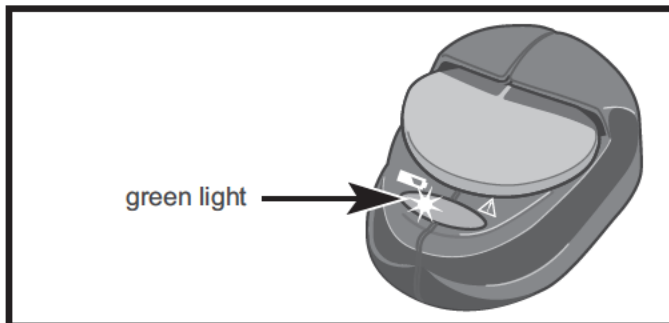

- 4 When charging is complete, the green light on the charger will stay on, without flashing, for 15 to 20 seconds and then turn off.
- 5 After the green charger light turns off, disconnect the transmitter from the charger. The green light on the transmitter will flash for approximately five seconds and then turn off.
- 6 If the green light on the transmitter does not flash, reconnect it to the charger for at least one minute.
- 7 After removing the transmitter from the charger, wait at least one minute before connecting it to a sensor or tester.

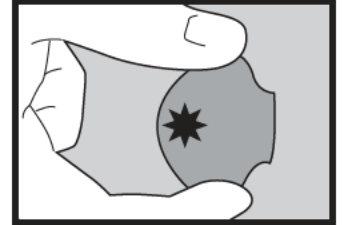

Inserting the sensor

This section shows you how to insert the sensor into the fatty layer under your skin so that it can measure the glucose in your body. You can insert the sensor manually or use the Sen-serter insertion device that came with your system. It is highly recommended that you use the Sen-serter to make sure that the sensor is placed correctly.

The sensor determines how much glucose is in the fatty layer under the skin at that site (location), then sends a signal with this information to the transmitter. The monitor translates the signal and displays a sensor reading on your monitor.

WARNING: The sensor may affect your medical conditions or medications. Bleeding, swelling, irritation and/or infection at the insertion site are possible risks associated with inserting the sensor and sometimes result from improper insertion and maintenance of insertion site. Please discuss these conditions, medications and/or reactions with your healthcare professional before using the sensor.

Before you begin

This section contains important information to prepare your sensor for insertion, where you should and should not insert the sensor, and how to prepare the insertion site.

Make sure that the transmitter is fully charged

Be sure you have set up and charged the transmitter before inserting the sensor.


Make sure that the sensor is ready to be inserted

If the sensor has been refrigerated, remove the sensor package from the refrigerator. To avoid condensation, make sure that you allow the sensor package to reach room temperature before opening it.

Where to insert the sensor

CAUTION: Never insert the sensor within 5.0 cm (2 inches) of an insulin pump infusion site, or within 7.5 cm (3 inches) of a manual insulin injection site.

Choose a site with an adequate fatty layer for sensor insertion. Be sure to rotate the sensor sites so that they do not become overused. Shown here are the best body areas (shaded) for sensor insertion. The clinical trials of sensor accuracy have been based on sensors inserted in the abdominal area. Sensor performance may differ when other insertion sites are used.

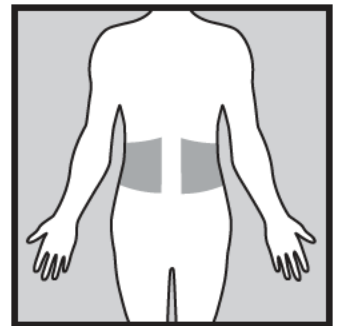

Do NOT insert the sensor in




- Frequently used injection or infusion sites.
- The 5.0 cm (2 inches) area around your navel.
- Sites where clothing rubs your skin or limits movement.
- Scarred or damaged tissue.
- Body sites that you move frequently.
- Sites on your body without enough fatty tissue.

How to prepare the insertion site

- Make sure that you clean the insertion site with alcohol and that the site is dry before you insert the sensor.
- Do NOT use sticky skin preparation solutions before you insert the sensor.
- You may use a sticky intravenous (I.V.) preparation solution after you insert the sensor, and before you apply a sterile dressing, to help the dressing stick to your skin.




Always refer to the instructions that came with your glucose sensor.

- 1 Wash your hands.
- 2 Clean the sensor site with alcohol. Let it air dry.
- 3 Remove the sensor from the package by holding the base or tape. Do not hold the sensor by the introducer needle handle.
- 4 Place the sensor in the Sen-serter until it fits snugly.
- 5 Holding the white tape as shown, remove the clear tape using a counterclockwise motion.
- 6 Place your fingers on the back of the white tape and push the carrier down until it clicks into place.

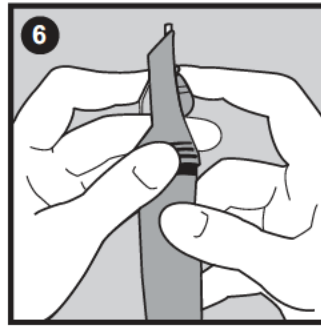
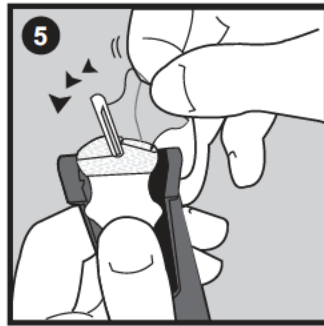
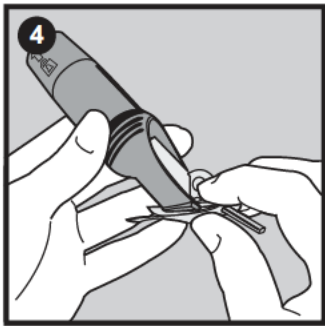


- 7 Turn the white button to lock the Sen-serter. Remove the needle guard from the introducer needle.
- 8 Rest the Sen-serter legs flat on your skin so the Sen-serter is at a *45-degree angle or greater* to the insertion site. Bleeding can occur if you insert the sensor at an angle less than 45 degrees. Place two fingers of your other hand on the Sen-serter legs to maintain the correct angle.
- 9 Turn the white button to unlock the Sen-serter. Press the white button to insert the sensor.

- 10 Make sure the sensor is inserted. If the sensor is not fully inserted, manually push it into place. The sensor should lay flat on your skin.

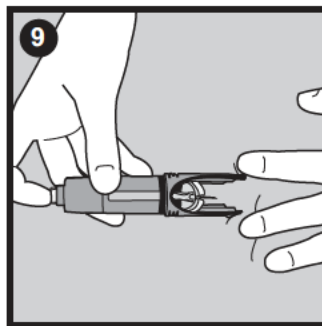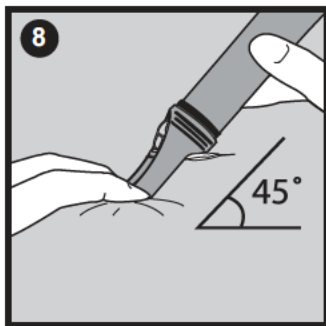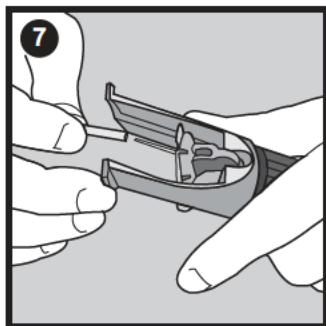

- 11 While holding the sensor in place, gently slide the Sen-serter away from the sensor. Do not twist, bend or lift the Sen-serter while removing it from the sensor.
- 12 While still holding the sensor in place, remove the white paper from the adhesive pad. Press the adhesive against your skin.
- 13 Hold the sensor with two fingers on the base, and gently remove the introducer needle. Do NOT rotate the introducer needle when removing. Dispose of the needle in sharps container.

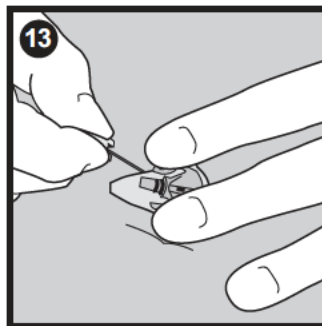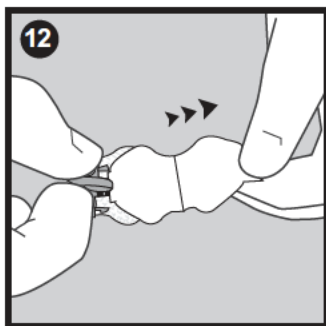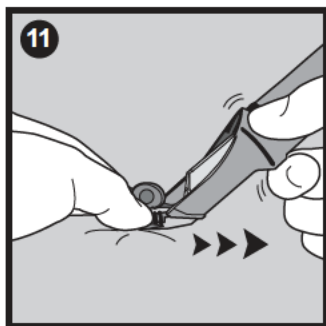

- 14 Wait five minutes after insertion before connecting the transmitter to the sensor to make sure the sensor is ready to communicate. Make sure that the site is not bleeding before connection. If bleeding occurs, apply pressure using a sterile gauze or clean cloth for three minutes.
- a. When bleeding stops, attach the transmitter to the sensor.


CAUTION: If bleeding does NOT stop, do NOT connect the transmitter to the sensor.

- b. If bleeding does not stop, do the following:
 - i Remove the sensor and discard.
 - ii Reapply pressure using a sterile gauze or a clean cloth until the bleeding stops.
 - iii Insert a new sensor in a different location.

Connecting the transmitter to the sensor

The transmitter must be set up and fully charged, and the sensor inserted, before you connect it to the sensor. If you connect the transmitter to a sensor that has not been inserted into your body, the transmitter will not flash or send signals to a monitoring device, and the transmitter battery will lose energy. After you disconnect the transmitter from the charger, allow at least one minute before you connect it to the sensor.

- 1 After the sensor is inserted, wait five minutes before connecting to the transmitter. Check for bleeding as described above, then connect the transmitter to the sensor.
- 2 Touch the inserted sensor at the back of the assembly to prevent movement.

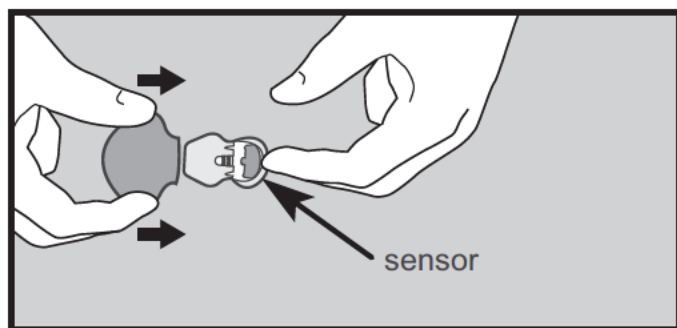

- 3 Hold the transmitter as shown and line up the two notches on it with the flexible side arms of the sensor. The transmitter's flat side should face the skin.
- 4 Push the transmitter onto the sensor arms until they snap into the transmitter notches. If properly connected, within 20 seconds the transmitter light will flash green for approximately 10 seconds.
- 5 If the transmitter light does not flash, disconnect it from the sensor, wait for one minute and then reconnect. If the transmitter light still does not flash, charge the transmitter, then repeat this procedure.

- 6 If the transmitter light flashes green after it is connected or reconnected to the sensor, then use your monitor to start communicating with the sensor and to start the sensor initialization. (See the next section for starting your sensor.)
- 7 After the transmitter successfully sends signals to the monitor, you may choose to put an occlusive dressing over the transmitter and the sensor.

Starting the sensor

This section shows you how to start the sensor so that it is ready to take glucose measurements and send them to the transmitter. You use the monitor to start the sensor.

To start the sensor:

- 1 Open the SENSOR START MENU.

MAIN MENU > Sensor > Sensor Start

The SENSOR START MENU shows **New Sensor**.

- 2 Press **ACT**. The NEW SENSOR screen shows a message reminding you to connect the transmitter to the sensor.
- 3 If you have not connected the transmitter to the sensor, connect it now (see previous section for the steps to use). However, if the NEW SENSOR screen times out and the HOME screen appears, DO NOT connect the sensor. Start this procedure again from the beginning.
- 4 Press **ACT**. The **SENSOR READY 2 HOURS** screen appears.
- 5 Press any button. The sensor is started.

After two hours, your monitor will prompt you to enter the first meter BG to calibrate the sensor.

Entering the first meter blood glucose

You cannot enter a meter blood glucose (meter BG) measurement into your monitor until you receive a METER BG NOW alert. If you enter a meter BG measurement before you receive this alert, it will not be used by your system for a calibration.

Your monitor will show a METER BG NOW alert after the two-hour initialization period. This alert prompts you to enter a meter BG reading into the monitor to calibrate the sensor. This calibration will be successful only if the blood glucose entry is in the range of 40-400 mg/dL (2.2 to 22 mmol/L).

To enter the first meter blood glucose measurement:

- 1 Press ESC, then ACT to clear the METER BG NOW alarm.
- 2 Check the sensor icon on the screen to make sure that the monitor and the transmitter are communicating properly.
 -  Communicating. Go to the next step and enter your first meter BG measurement.
 -  Not communicating. Do not enter your first meter BG measurement. Check the following:
 - The transmitter is connected correctly to the sensor.
 - The sensor has been started.
 - The monitor is not too far from the transmitter.
 - If there is still no communication, see the *Troubleshooting and alarms/alerts* chapter for information on how to restore the communication between the monitor and the transmitter.
- 3 Take a meter BG measurement. **Immediately enter the meter BG measurement into your monitor by following the rest of the steps in this procedure.**
- 4 Open the ENTER BG screen.
MAIN MENU > Sensor > Enter Meter BG

The ENTER BG screen shows three dashes.
- 5 Enter the number from your meter (the value on the screen must match your meter measurement). Values below 40 mg/dL (2.2 mmol/L) or above 400 mg/dL (22.2 mmol/L) cannot be used for calibration.
- 6 Press ACT. The system is calibrated.

After you enter your first meter BG measurement, you are ready to use your system. You must calibrate your system regularly.

Meter option

OneTouch® UltraLink™ and Paradigm Link™ glucose meters can be programmed to communicate wirelessly to the Guardian® REAL-Time System. When the wireless feature is turned on, glucose values transmitted from these meters to the Guardian monitor are stored in memory and used for glucose sensor calibration.

Your monitor is set at the factory with the meter option turned off. If you have turned your meter option on, you should turn off the wireless feature when using your glucose sensor. By doing so, you can manually enter calibration glucose values when glucose is stable and not changing rapidly.

Instructions for turning the wireless feature **OFF** are outlined below.

- From the main menu, select **UTILITIES** and press **ACT**
- Scroll to **METER OPTION** and press **ACT**
- Select **OFF** and press **ACT**

You can set up your monitor to wirelessly receive your blood glucose reading from the OneTouch® UltraLink™ Meter or the BD Paradigm Link® glucose meters (see the *Meter rules* section in this chapter). These meters may not be available in all countries. Check with your local Medtronic Diabetes representative. Programming a meter ID into your monitor links your monitor to the meter and allows fingerstick blood glucose calibration values to be automatically transmitted to your system. If you do not link a meter to your monitor, you must enter your blood glucose readings manually. Each meter has its own unique ID. You can link up to three meters to your monitor. The use of RF meters with the monitor reduces monitor battery life.

When the monitor is at the HOME screen, it will beep or vibrate when it receives a blood glucose measurement from the meter. The measurement will appear on the screen.

You have to turn on the meter option to add, delete or review the meter ID(s) programmed in your monitor. The meter ID is the serial number printed on the back of the meter. See the user guide that came with your meter for detailed information on how to use it.

Meter rules

If you want your monitor to communicate with the meter, the following conditions must apply:

- 1 The meter option must be turned on and programmed. Refer to the instructions in this section.
- 2 Your monitor must be within 1.2 meters (4 feet) of your meter to receive the blood glucose reading.
- 3 The monitor cannot have a LOW BATTERY alert condition.
- 4 The blood glucose measurement from the meter will appear as the default blood glucose value on the ENTER BG screen. The monitor will not display a reading that is older than 12 minutes on the ENTER BG screen.
- 5 Do not use the RF meter to send your blood glucose readings to the monitor while onboard aircraft. Manually enter your blood glucose.

Add, delete, review meter IDs

Select **Meter Options** in the Utilities Menu screen when programming your meter.

If you are not sure that your meter ID is entered in your monitor check the REVIEW METER ID screen.

You have to turn on the meter option to add, delete, or review the meter ID(s) programmed in your monitor

- 1 Go to the METER OPTION screen.

Main > Utilities > Meter Options > Meter Option

- 2 Select **On**, then press **ACT**. The METER ID MENU appears.
- 3 Add, delete or review your meter ID(s) as desired.

Adding meter IDs

- a. Select **Add ID**, then press **ACT**.
- b. Use the up and down arrow buttons to enter each of the six ID characters. Press **ACT** after each entry.
- c. After you set the last character of the ID, the screen will return to the METER ID MENU.

Deleting meter IDs

- a. Select **Delete ID**, then press **ACT**.
- b. Select the meter ID that you want to delete and press **ACT**.
- c. The selected ID is now deleted.

Reviewing meter IDs

- a. Select **Review ID**, then press **ACT**.
 - b. The programmed IDs will show in the REVIEW METER ID screen.
- 4 Exit the menus when you are done.

Using your system

In this chapter:

- Calibrating your system
- Capturing events
- Reading the sensor glucose graphs
- How to check for rapid changes in sensor glucose
- Checking the status of the system
- How to make glucose alerts silent
- Viewing the sensor update history
- Viewing your current settings
- Disconnecting the transmitter and removing the sensor
- Using your system in water

Calibrating your system

Calibrating your system helps to make sure that the glucose readings from the sensor are as accurate and reliable as possible. You must periodically enter meter BG measurements to calibrate your system; however, do not enter a meter BG until at least two hours after you have started a new sensor.

CAUTION: Sometimes your glucose levels change so rapidly that blood glucose measurements taken during those times should not be used to calibrate your system. It is recommended that you do not calibrate your system when your blood glucose may be changing rapidly, such as: within the first hour after eating; when rapid change arrows appear on your monitor; or if you recently received a Rate of Change alert. It is recommended that you calibrate your system only when your blood glucose is stable, such as before meals or going to bed.

You may negatively affect your calibration accuracy if you allow your meter to automatically transmit fingerstick blood glucose calibration values to your system when your blood glucose may be changing rapidly. If you think your blood glucose may be changing rapidly, check the sensor glucose (SG) information on your monitor before taking a fingerstick BG measurement. If your SG shows rapid changes, unlink your meter before taking your fingerstick BG measurement.

Guidelines for calibrating the system

This section includes important recommendations and information about calibrating the system.

When should I calibrate the system?

When you need to calibrate your system, the monitor will show a METER BG NOW alert at the following times:

- At least twice a day (once every 12 hours), whether or not you have inserted a sensor that day.
- Two hours after you insert a new sensor.
- Within six hours of first calibrating a new sensor, a failed calibration, or an invalidated calibration (for example, entering a meter BG for calibration, then removing the battery for more than ten minutes).

What are the best times to calibrate the system?

It is recommended that you calibrate your system when your blood glucose is **not** changing rapidly, such as:

- Before meals
- Within the first hour after eating a meal or snack
- Before going to bed

What should I do if I have my RF meter is set to automatically transmit my fingerstick BG measurement and my blood glucose is changing rapidly when I take the measurement?

It is recommended that you check the sensor glucose (SG) information on your monitor before taking a fingerstick BG measurement. If the SG information shows rapid changes, then turn off the wireless communication feature before taking your fingerstick BG measurement. After you take the BG measurement, you may turn the wireless communication feature back on.

What happens if I do not calibrate often enough?

This example assumes that you are using the Cal Reminder default setting of one hour.

Time since last calibration	What happens
At least 11 hours	If you have not changed the default Cal Reminder setting (one hour), the system reminds you to calibrate the system by sending a METER BG BY alert one hour before your next calibration is due. You can change the CAL REMINDER alert setting so that you receive an alert earlier than the default.
At least 12 hours	• The system sends you a METER BG NOW alert. • The system stops taking real-time glucose measurements. • The system stops sending alerts when you reach your High/Low Glucose limits. • The system will not take glucose measurements again until you enter a new Meter BG measurement.

To calibrate your system:

- 1 Check the sensor icon on the screen to make sure monitor and transmitter are communicating.
 -  Not communicating. Do not enter your meter BG measurement for calibration. Check the following:
 - The sensor and transmitter are correctly connected.
 - The sensor has been started.
 - The monitor is not too far from the transmitter.
 - The transmitter is charged.
 - The correct transmitter ID has been entered.
 -  Communicating. Go to the next step to calibrate your system.

If there is still no communication, see the *Troubleshooting and alarms/alerts* chapter for information on how to restore communication between the monitor and the transmitter.

- 2 Take a meter BG measurement. **Immediately enter the meter BG measurement into your monitor by following the rest of the steps in this procedure.**
- 3 Open the ENTER BG screen.

Main > Sensor > Enter Meter BG

The ENTER BG screen shows three dashes.

- 4 Select the number from your meter. The number on the screen must match the measurement from your meter, and values below 40 mg/dL (2.2 mmol/L) or above 400 mg/dL (22.2 mmol/L) cannot be used for calibration.
- 5 Press **ACT**. The calibration has been successfully entered. It takes 10-15 minutes after entering your meter blood glucose value for the measurement to be accepted, and for the system to show those values.

If you enter an incorrect calibration value:

- 1 If you enter an incorrect calibration value, immediately re-enter the correct value (if two values are entered within 15 minutes of one another, the system uses the most recently entered value).

Capturing Events

The Capture Event feature allows you to electronically save certain types of information that is used with Medtronic CareLink Therapy Management software (not available in all locations) to generate treatment reports that you can share with your healthcare professional. Some examples include:

- Non-calibration meter BG measurements (used to manage your diabetes but not to calibrate the system)
- The type and amount of insulin you use
- The amount of carbohydrates you consume
- The intensity and duration of exercise you do
- Other (this can be used for anything else as needed)

Make sure you enter events when they happen. The system records the time of the entry and cannot be changed to show an earlier time. You also cannot change entries after you have put the information into your monitor (except meter BG entries used to calibrate the sensor).

Capture Event also allows you to manually enter and save the following information, which can also be sent to Medtronic CareLink Therapy Management software (not available in all locations):

- Sensor glucose measurements
- Sensor glucose alert settings
- Sensor glucose limit alerts
- Meter blood glucose measurements that are *used to calibrate the system*

Entering non-calibration meter BG measurements

Non-calibration meter BG measurements are those taken when you will not use them for calibration, such as when you eat or when your blood glucose is rising or falling rapidly.

- 1 Make a note of the non-calibration meter blood glucose measurement.
- 2 Open the ENTER BG screen.

Main > Capture Event > Enter BG

The ENTER BG screen flashes with the most recent meter BG measurement you entered.

- 3 Enter your new measurement, then press **ACT**. A message asks if you want to save the new measurement that is shown on the screen. The **Yes** option is selected.
- 4 Make sure that the number shown on the screen is correct.
 - If the information is correct, press **ACT**.
 - If the information is not correct, select **No**, then press **ACT**. The CAPTURE EVENT menu reappears with **Enter BG** selected. Repeat the procedure and enter the correct information.
- 5 The BG TO UPDATE SENSOR? screen appears. This allows you to use this meter BG reading to calibrate your sensor if you choose. For this screen to appear, the sensor must be turned on, started and initialized, and communicating with the transmitter. Your meter BG measurement must be between 70 mg/dL (3.8 mmol/L) and 250 mg/dL (13.8 mmol/L).
 - To update the sensor with this meter BG reading, select **Yes**.
 - If you do not want to update the sensor with this reading, select **No**.
 - For the Update Sensor screen to appear as described above, the sensor must be turned on, started and initialized, and communicating with the transmitter. Your BG measurement must be between 70 mg/dL (3.8 mmol/L) and 250 mg/dL (13.8 mmol/L).

Entering insulin injection information

Follow these steps to enter the type and amount of insulin you use.

- 1 Make a note of the amount and type of injection you are taking.
- 2 Open the INSULIN TYPE screen.

Main > Capture Event > Insulin Marker

The INSULIN TYPE screen automatically selects the most recent insulin type you entered.

- 3 Select the correct insulin type for the injection. The insulin types are:
 - Fast Acting
 - Short Acting

- Intermediate Acting
 - Long Acting
 - Mixed 70/30
 - Mixed 50/50
- 4 Press **ACT**. The ENTER INSULIN screen flashes the insulin amount that you most recently entered.
 - 5 Enter the amount of the injection, then press **ACT**. A message appears asking you if you want to save the new information that is displayed on the screen. The **Yes** option is selected.
 - 6 Make sure that the type and amount of insulin shown on the screen is correct.
 - If the information is correct, press **ACT**. The information you entered is saved to the system and can now be used in reports.
 - If the information is **not** correct, select **No**, then press **ACT**. The CAPTURE EVENT menu shows with **Enter BG** selected. Repeat the procedure to enter the correct information.

Entering information about carbohydrates

This section shows you how to enter information about the carbohydrates you eat or drink with meals or snacks.

Before you begin

Before you can enter carbohydrate information, you must select the carbohydrate units (Carb Units) for the Capture Event feature. You can set the carbohydrate units to grams or exchanges.

To set the carbohydrate units:

- 1 Open the CARB UNITS screen.
Main > Capture Event > Carb Units
 The CARB UNITS screen shows **Grams** selected.
- 2 Select either **Grams** or **Exchanges**, then press **ACT**. The CAPTURE EVENT menu shows the carbohydrate unit you selected.

To enter carbohydrate information:

- 1 Determine the total units of carbohydrates in the meal, snack or drink that you plan to consume.
- 2 Open the ENTER MEAL screen.
Main > Capture Event > Meal Marker
 The ENTER MEAL screen shows **Breakfast** selected.
- 3 Select the correct meal type for what you consume: Breakfast, lunch, dinner or snack.

- 4 Press **ACT**. The ENTER FOOD screen flashes either with three dashes, or with the number of carbohydrate units you last entered.
- 5 Enter the carbohydrate units, then press **ACT**. A message asks if you want to save the information on the screen. The **Yes** option is selected.
- 6 If the number on the screen is correct, press **ACT**. The information you entered is saved and can be used in reports. If the information is **not** correct, select **No**, then press **ACT**. The CAPTURE EVENT menu shows with **Enter BG** selected. Repeat the procedure and enter the correct information.

Entering exercise information

Follow these steps to enter information about the exercise you do.

- 1 Make a note about how long you exercised (duration), and how difficult or easy the exercise was (intensity).
- 2 Open the EXERCISE INTENSITY screen.

Main > Capture Event > Exercise Marker

The EXERCISE INTENSITY screen flashes the intensity type you selected during your last entry.

- 3 Select the intensity type. The intensity types are:
 - High (difficult)
 - Medium (some difficulty)
 - Low (easy)
- 4 Press **ACT**. The EXERCISE DURATION screen flashes the length of time you entered for your most recent exercise marker.
- 5 Enter the length of time you exercised, then press **ACT**. A message asks if you want to save the information that is displayed on the screen. The **Yes** option is selected.
- 6 If the number on the screen is correct, press **ACT**. The information you entered is saved and can be used in reports. If the information is **not** correct, select **No**, then press **ACT**. The CAPTURE EVENT menu shows with **Enter BG** selected. Repeat the procedure and enter the correct information.

Entering Other markers

This section shows you how to enter markers other than Enter BG, Insulin Marker, Meal Marker, or Exercise Marker. When you enter Other markers, you only enter the type of marker (Other); the time is automatically recorded by the system. Examples of Other markers include: when you take medications, when you feel ill, when you are under stress, etc. These markers then create reports that appear with the name of the marker (Other) and the time you entered it.

To enter Other markers:

- 1 Go to the OTHER screen.

Main > Capture Event > Other

A message appears: Do you want to enter an Other marker?

- 2 Select Yes, then press ACT. Your monitor saves your Other marker and returns to the CAPTURE EVENT screen.

Reading the sensor glucose graphs

Your monitor shows an updated, real-time glucose measurement every five minutes. This measurement is generated by data sent from the sensor to the transmitter, and then to the monitor, where it is converted to a glucose measurement shown as 3-, 6-, 12- and 24-hour graphs. Make sure your graph screen does not say DEMO in the upper left corner. If it does, you must turn the Sensor Demo off. The graphs include the following information:

- The most recent sensor glucose or the reason that a sensor glucose measurement is not displayed.
- The historical sensor glucose measurements taken for the last 3, 6, 12 or 24 hours, including the value and the time the sensor glucose measurement was taken, or the reason it is not displayed.
- Arrows that show the relative rate at which the most recent sensor glucose levels are rising or falling.

Your system will not alert every time an SG measurement is not taken, so your alarm/alert history may not match the number of times your graphs do not show an SG measurement. The following conditions can prevent a real-time sensor glucose measurement from being taken:

- A Lost Sensor alert.
- A Sensor Error alert.
- A Bad Sensor alert.
- A new sensor that you just inserted is being initialized and has not yet started taking sensor glucose measurements.
- A Cal Error alert.
- More than 12 hours has passed since the last system calibration.

- Find Lost Sensor (15 minutes).
- Weak Signal.
- Reconnect old sensor (two hours).

Two other conditions also prevent a real-time SG measurement and generate a WARM UP message on your sensor glucose graphs:

- A new sensor that you just inserted is being initialized and has not yet started taking SG measurements.
- More than 12 hours have passed since the last calibration.


Opening and viewing the graphs

The following table describes how to use the graphs.


If you want to:	Do the following:
Open the graphs	From the HOME screen, press ESC . The last graph you viewed (3-, 6-, 12- or 24-hour) shows the details for your most recent sensor glucose (SG) measurement.
View information in a graph	If there is no SG measurement at the time you select, the reason for the missing measurement will show on the right side of the graph (for example, Lost Sensor). • When you open the graphs, the cursor (flashing vertical line) is always on the right edge of the graph and the graph shows the most recent SG measurement. • Press  to move the cursor left to see an earlier SG measurement. • Press  to move the cursor right to return to more recent SG measurements.
View different graphs	• With the cursor at the far right edge of a graph, you can press  to open the next graph. • If the cursor is anywhere to the left of the far right edge, press ESC to move it to the right, then press  to open the next graph.

The graphs

All graphs show High and Low Glucose Limits lines, a real-time SG measurement line, a data section, and the cursor (flashing vertical line). When you open any graph, the cursor flashes on the right edge of the graph, the data section shows the most recent SG measurement (or the reason why no measurement shows), and the time in the data section matches the time at the top of the screen.

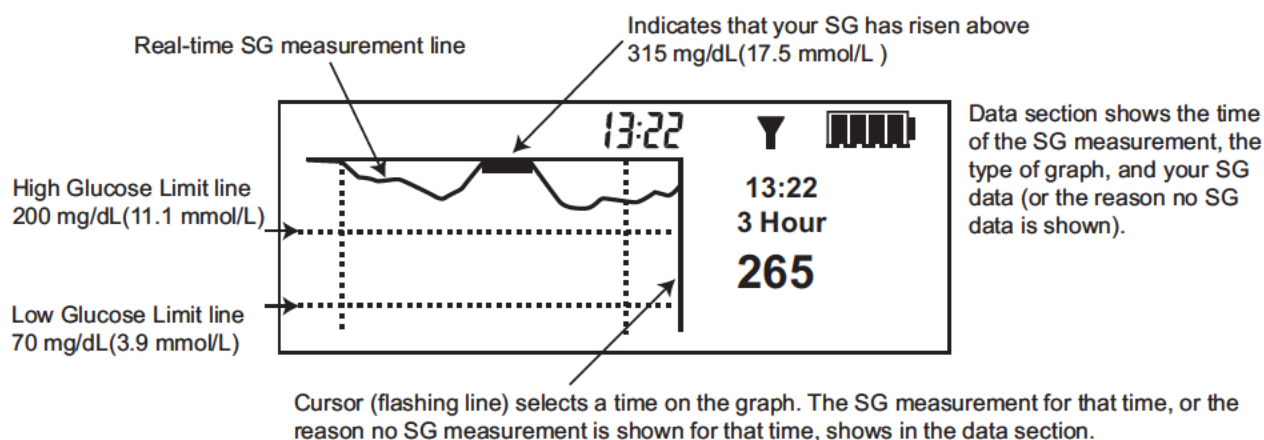

When you move the cursor left to select an earlier SG measurement, the data section of the graph turns black, the word History appears, and the time in the data section changes to show the time when the SG measurement was taken (or the reason why no SG measurement shows for that time). Because your system does not alert every time there is no SG measurement, your alarm/alert history may not match the number of times your graphs have SG measurement gaps.

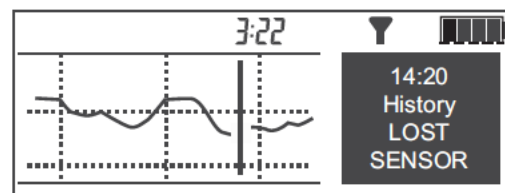


Examples of real-time sensor glucose graphs

The following pages show examples of the different graph types.


3-hour graph

Each part of the graph between the vertical dashed lines represents one hour. In this example, the selected SG measurement was taken at 14:50 (2:50 p.m.), and the measurement at that time was 150 mg/dL (8.3 mmol/L).

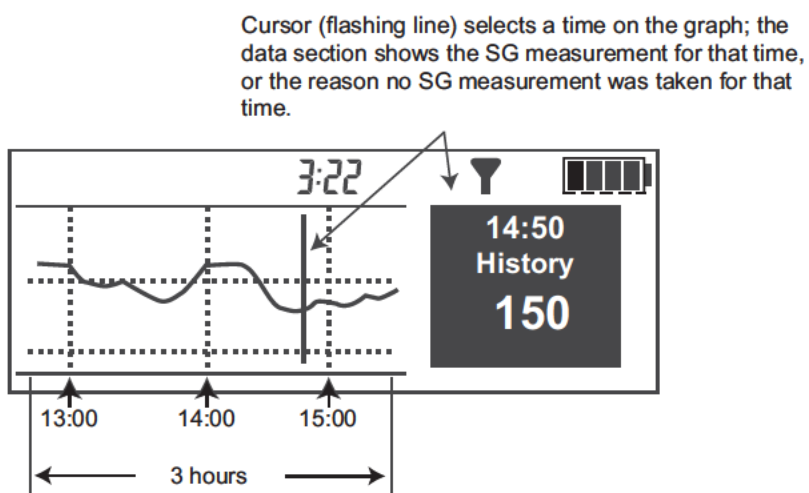


6-hour graph

Each part of the graph between the vertical dashed lines represents one hour. In this example, the selected SG measurement was taken at 1:12 a.m., and the measurement at that time was 149 mg/dL (8.2 mmol/L).

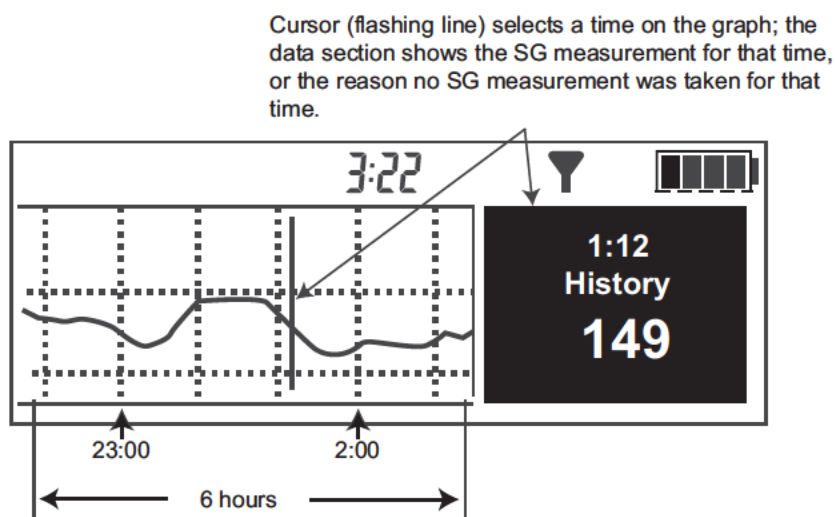

12-hour graph

Each part of the graph between the vertical dashed lines represents one hour. In this example, LOST SENSOR is the reason no measurement is shown, and it occurred at 14:20 (2:20 p.m.).

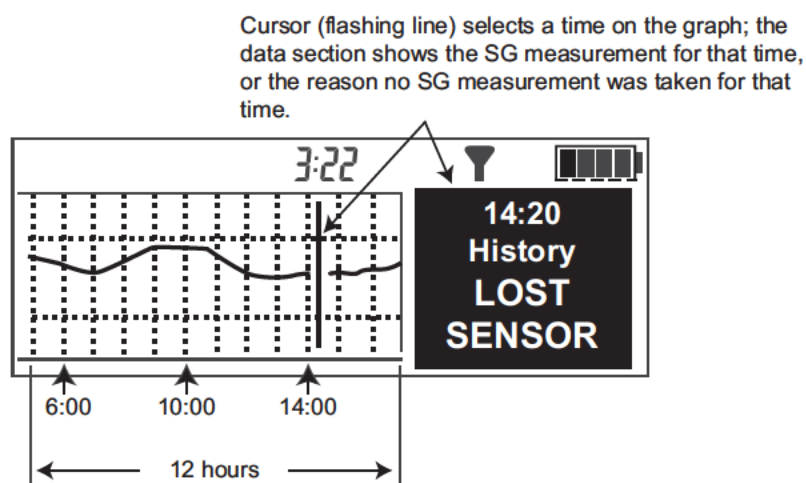


24-hour graph

Each part of the graph between the vertical dashed lines represents six hours. The shaded area helps you more easily see the sensor glucose measurements for the previous night. In this example, the selected SG measurement was taken at 11:27 a.m. and the measurement at that time was 119 mg/dL (6.6 mmol/L).

The shaded area is the 12 hours from 18:00 to 6:00.

Cursor (flashing line) selects a time on the graph; the data section shows the SG measurement for that time, or the reason no SG measurement was taken for that time.

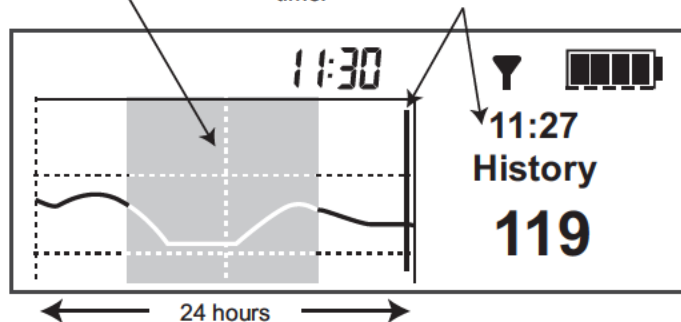


How to check for rapid changes in sensor glucose

If your sensor glucose (SG) rises or falls faster than a certain rate, Rapid Change arrows automatically appear in the graphs next to the SG number, as shown in this section. You can also set your Rate of Change alerts to tell you when these arrows appear on your SG graphs.


The rapid change arrows

These examples show how the Rapid Change arrows tell you when your SG rises or falls faster than a certain per-minute rate. The arrow direction shows if your SG measurements are rising or falling. The number of arrows (one or two) shows how rapidly your SG is changing.


- One up arrow ↑ shows that your SG has been rising at a rate of 1.0 to 2.0 mg/dL (0.06 and 0.11 mmol/L) per minute.

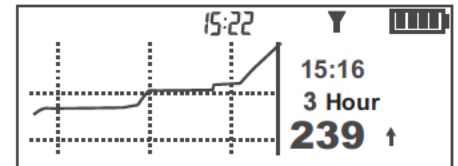


- One down arrow ↓ shows that your SG has been falling at a rate of 1.0 to 2.0 mg/dL (0.06 and 0.11 mmol/L) per minute.

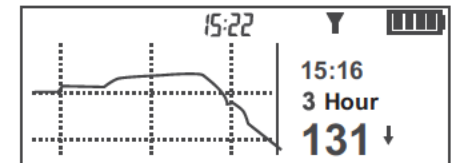

- Two up arrows $\uparrow\uparrow$ show that your SG has been rising at a rate of 2.0 mg/dL (0.11 mmol/L) or more per minute.

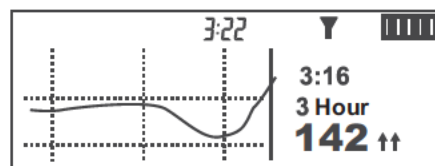

- Two down arrows $\downarrow\downarrow$ show that your SG has been falling at a rate of 2.0 mg/dL (0.11 mmol/L) or more per minute.

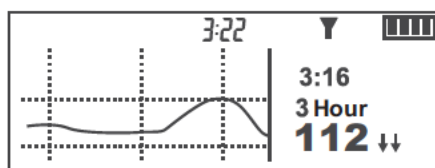

Checking the status of the system

The system allows you to check its status (condition) in the following areas:

- Condition of the batteries.
- The age of the sensor.
- Connection and communication between the sensor and transmitter, and the transmitter and monitor.
- When the next meter BG measurement is needed to calibrate the system.

Checking the monitor battery status

The battery status shows you how much energy remains in the battery.

- 1 Make sure the HOME screen appears. If it does not, press **ESC** until it appears.
- 2 Press the **ESC** button until the STATUS screen appears.

One of the following status types is displayed:

Status	Description
Normal	The battery has enough energy to run the features of the system.
Low	The battery is running low on energy. The battery should be replaced soon.
Off	The battery is out of energy. The system shuts down automatically. The battery should be replaced as soon as possible.

Checking the sensor status

The sensor status tells you when your next calibration is due, and how long the sensor has been used.

Before you begin

The sensor and transmitter must be connected, and the sensor must be turned on, before you can check the sensor status.

To check the sensor status:

- 1 From the HOME screen, press **ESC** until the SENSOR STATUS screen shows the following status information:

Status	Description
Next Cal	Shows when your next calibration is due.
Sensor Age	The total number of days and hours that the sensor has been used.
Sensor Isig	This information is about the signal of the sensor, not sensor status, and is used only by the Medtronic HelpLine for diagnostic purposes.

Checking the transmitter status

The transmitter status tells you if the transmitter is connected to the sensor, and the status (condition) of the transmitter battery.

To check the transmitter status:

- 1 From the HOME screen, press **ESC** until the SENSOR STATUS screen shows.
- 2 Press  to see the following transmitter status information:

Status	Description	
Transmtr Batt	Shows three possible battery conditions:	
	Good	The transmitter battery has enough energy to run all features of the transmitter.
	Low	The transmitter battery is running low on energy. The transmitter should be re-charged within 24 hours. It will continue to send data until the battery is completely out of energy.
	Bad	The transmitter battery is out of energy. The transmitter shuts down automatically, and no sensor data is sent to the monitor. The transmitter should be recharged as soon as possible.
Transmtr Ver	The version number of the transmitter.	
Transmtr ID	The serial number of the transmitter.	

How to make glucose alerts silent

The Alert Silence feature allows you to make glucose alerts silent for a pre-set amount of time between 30 minutes and 24 hours.

You do not lose information when you silence the glucose alerts. The system still records each glucose alert time and SG measurement.

These are your Alert Silence options:

- **Off** — This means the alerts are **On**: the monitor will beep or vibrate when any sensor glucose alerts occurs.
- **Hi only** — The monitor will not beep or vibrate when the High Glucose Limit alert occurs during the specified time.
- **Lo only** — The monitor will not beep or vibrate if the Low Glucose Limit alert occurs during the specified time.
- **Hi and Lo** — The monitor will not beep or vibrate if the High Glucose Limit or Low Glucose Limit alert occurs during the specified time.
- **All** — The monitor will not beep or vibrate if any sensor glucose alert occurs during the specified time.

To make alerts silent:

- 1 Open the SENSOR MENU.

Main Menu > Sensor

- 2 The SENSOR MENU shows you the current status of the Alert Silence feature.
- 3 Depending on your current settings, one of the following **Alert Silence** options will appear in the SENSOR MENU under the **Sensor Update Hist.** option.
 - Alert Silence: Off
 - Hi only: Silent
 - Lo only: Silent
 - Hi/Lo: Silent
 - All: Silent (Includes the Low Battery alert for the monitor battery)
- 4 To change the current alert silence settings, select **Alert Silence**, then press **ACT**.
The ALERT SILENCE screen shows the alert (or alerts) that are set to be silent already selected.
- 5 Select the alert silence option you want, then press **ACT**. The SET DURATION screen flashes one of the following:
 - A set of dashes (if this is the first time you have used the Alert Silence feature, or you just cleared your settings).
 - The length of time of the last alert silence that you entered.

- 6 Select your alert silence length between 30 minutes (minimum) and 24 hours (maximum), then press **ACT**.
The **SENSOR MENU** shows the new alert that you selected to be silent.
- 7 Press the **ESC** button until the **MAIN MENU** appears.

Viewing the sensor update history

A sensor update is a meter BG measurement that is entered to calibrate the system.

Your system automatically records manually entered meter BG measurements. You can review the last **28** sensor updates. The **date**, **time** and **meter blood glucose measurement** is recorded for every meter BG measurement that is successfully used to calibrate the system. It takes approximately 20 minutes from the time your meter BG measurement is entered and accepted until it shows in the **SENSOR UPDATE HISTORY** screen.

To review the Sensor Update History:

- 1 Open the **SENSOR UPDATE HISTORY** screen.

Main Menu > Sensor > Sensor Update History

The **SENSOR UPDATE HISTORY** screen shows the recent meter BG measurements that were used to calibrate the system.

- 2 The date and time of the calibration occurs 15 minutes after the meter BG measurement is entered.
- 3 Use the  and  buttons to view the information on the screen. Your most recent sensor update is displayed at the top of the list. If a calibration was not successful, it will not display in the **SENSOR UPDATE HISTORY** screen.

Viewing your current settings

You can review your current settings at any time. This is useful for making sure that you use the correct settings when troubleshooting and for checking a setting before changing it.

- 1 From the HOME screen, open the REVIEW SETTINGS screen.

Main Menu > Sensor > Sensor Setup > Review Settings


The REVIEW SETTINGS screen appears.

- 2 Press the  and  buttons to scroll through your settings.
- 3 The High Snooze setting will not appear unless you switch on at least one of the following alerts - High Glucose, Rise Rate, or High Predictive. The Low Snooze will not appear unless you switch on one of the following alerts: Low Glucose, Fall Rate, or Low Predictive.

Disconnecting the transmitter and removing the sensor

Disconnecting the sensor from the transmitter

If you are not going to replace the sensor, turn the Sensor features **Off** to avoid getting a LOST SENSOR alert.

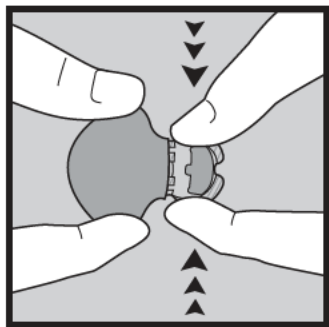

- 1 Hold the transmitter as shown, and pinch the flexible side arms of the sensor between your thumb and forefinger.
- 2 Gently pull the transmitter away from the sensor assembly. If you need to connect the transmitter to the tester, wait one minute after disconnecting the transmitter from the sensor.

Removing the sensor

Gently pull the sensor from your body to remove it. Place it in a sharps container.

Using your system in water

Your monitor is NOT waterproof: remove it before doing any water activities.

After your transmitter and sensor are connected, they form a water-tight seal to a depth of 8 feet (2.4 meters) for up to 30 minutes. You can shower, bathe or swim for up to 30 minutes without damaging them. However, hot water may decrease sensor life.

74 Using your system

Utilities

In this chapter:

- The data management feature
- The Alarm Clock
- The Block feature
- The Lock Keypad feature
- Running the Selftest
- User Settings

The data management feature

The system keeps a daily record of the following types of information for up to 31 days:

- **Daily totals** – Daily totals of certain information tracked automatically for you by the system.
- **Averages** – Averages of certain information calculated automatically for you by the system. There are two types:
 - An average over a number of days (you select the number of days)
 - An average for a single day

Both the SG information and meter BG information include daily totals and averages, as described in the following sections.

Sensor glucose information

The following table lists the different types of SG information that the system automatically records. The SG information is calculated by using all of the SG measurements that were received by the monitor during the day.

Item	Description
SG Average	The average sensor glucose measurement.
Above Hi Limit	The percentage of sensor glucose measurements that were above your High Glucose Limit.
Within Limits	The percentage of sensor glucose measurements that were between your Low Glucose Limit and your High Glucose Limit.
Below Low Limit	The percentage of sensor glucose measurements that were below your Low Glucose Limit.
SG Std. Dev.	The standard deviation of the sensor glucose measurements.
AUC Hi/Lo.	The High and Low AUC reference values. This is a reference value that provides a relative indication of the overall extent and duration of the glucose excursions on a given day. These numbers are shown in mmol/L or mg/dL, depending on your BG Unit setting.
#Hi/Lo Alarm	The total number of glucose alerts. The number to the left is the total number of high glucose alerts. The number to the right is the total number of low glucose alerts.

Meter blood glucose information

The table below lists the different types of meter BG information that the system automatically records. The meter BG information is calculated by using all of the automatic and manual meter BG measurements received by the monitor, including measurements that were manually entered into the monitor. You may view the data for a single day, or for a selected number of days.

An asterisk (*) next to a date means that the data for that day will not be used for the calculations that are displayed in the Daily Average option. Anything, such as a time or date change, that shortens a day at least one hour and twelve minutes (01:12) will cause the asterisk to display.

Item	Description
BG Average	The average meter blood glucose measurement.
Number BGs	The total number of meter blood glucose measurements taken.
Meter L/H	The lowest and highest meter blood glucose measurements used to calibrate the system (L is the lowest, and H is the highest).
Manual L/H	The lowest and highest manually entered blood glucose measurements (L is the lowest, and H is the highest). These can be used for calibration or non-calibration purposes.

To view data for a single day:

- 1 Open the DAILY TOTALS screen.
MAIN MENU > Utilities > Daily Totals
- 2 The BG AVG (Blood Glucose Average) screen displays:
 - To view totals for today, select **Today**.
 - To view totals for a previous day, select the date.
- 3 Press **ACT**. The SENSOR STATS screen shows the sensor glucose information for the day you selected (the day you selected is shown on the screen).
- 4 Press  to view more information.
- 5 To view meter BG measurement information for the same day you have selected in step 2, press **ACT** from the SENSOR STATS screen. The SENSOR STATS screen shows the meter BG information for the day you selected (the day you selected is shown on the screen).

- 6 Press  to view all of the information in the SENSOR STATS or METER STATS screen.
- 7 Continue viewing daily totals by following the steps above, or press **ESC** until the HOME screen appears.
The **ACT** button allows you to toggle between the SENSOR STATS and METER STATS screens. Pressing **ESC** from the METER STATUS screen will take you to the SENSOR STATUS screen.

To view data over a number of days:

- 1 Open the DAILY TOTALS screen.
MAIN MENU > Utilities > Daily Totals
The DAILY TOTALS screen shows with **Daily Average** selected.
- 2 Press **ACT**. The DAYS TO AVERAGE screen flashes the number of days you selected the last time you used this feature.
- 3 Select the number of days that you want to average. For example, if you want to see the averages for the last seven days, select 7.
- 4 Press **ACT**. A message tells you that the system is calculating the averages.
Then the SENSOR STATS screen shows the SG information for the number of days you selected. The number of days is also shown on the screen.
- 5 Press  to view more information.
- 6 To view meter BG measurement information for the same days you selected in step 3, press **ACT**. The METER STATS screen shows the meter BG information.
- 7 Press  to view more information.

The Alarm Clock

The Alarm Clock allows you to set up daily Alarm Clock reminders that help you remember to check your blood glucose, eat, calibrate the system, or do any other activity that you use to manage your diabetes. The default setting for this feature is **Off**.

When you set up an Alarm Clock reminder, you set the time when you want the system to send you the reminder. You can set a total of eight daily reminders.

The Alarm Clock reminder beeps or vibrates just like an alarm or an alert, and ALARM CLOCK shows on the screen. After the reminder shows on your screen, press **ESC**, then press **ACT** to clear it.

To set up daily Alarm Clock reminders:

- 1 Open the ALARM OPTION screen.
MAIN MENU > Utilities > Alarm Clock
- 2 Select **On/Set**, then press **ACT**. The ALARM CLOCK screen appears.
- 3 Select **Add Alarm**, then press **ACT**. The ADD ALARM screen shows a default time of 00:00.
- 4 Set the time you want for the reminder.
- 5 When using the 12-hour time setting, check the ADD ALARM to make sure that you have also selected either A for A.M. or P for P.M.
- 6 Press **ACT**. The reminder is set.
- 7 Repeat these steps to set up additional Alarm Clock reminders.
- 8 Press **ESC** until the HOME screen appears.

Reviewing your current Alarm Clock reminders

This procedure shows you how to view your current Alarm Clock reminders. The dashes in this example show that you may set up additional reminders. The total number of reminders that you can set up is eight.

- 1 Open the ALARM OPTION screen.
MAIN MENU > Utilities > Alarm Clock
The ALARM OPTION screen appears.
- 2 Select **On/Set**, then press **ACT**. The ALARM CLOCK screen appears.
- 3 Select **Review Alarms**, then press **ACT**. The REVIEW ALARMS screen shows current Alarm Clock reminders.
- 4 Press **ESC** until the HOME screen appears.

Deleting Alarm Clock reminders

This procedure shows you how to delete Alarm Clock reminders.

- 1 Open the ALARM OPTION screen.
MAIN MENU > Utilities > Alarm Clock
The ALARM OPTION screen appears.
- 2 Select **On/Set**, then press **ACT**. The ALARM CLOCK screen appears.

- 3 Select **Delete Alarms**, then press **ACT**. The DELETE ALARMS screen shows the first Alarm Clock reminder selected (highlighted).
- 4 Select the reminder you want to delete, then press **ACT**. The reminder is deleted.
- 5 Repeat these steps to delete additional Alarm Clock reminders.
- 6 Press **ESC** until the HOME screen appears.

The Block feature

When this feature is turned on, anyone who may have access to your monitor cannot set or change any of your settings. The only features that can be used are the ones that you use to view information. The vibrate alert type is disabled when Block is on. Discuss what uses and settings are best for you with your healthcare professional.

To turn the on Block feature:

- 1 Open the BLOCK OPTION screen.
MAIN MENU > Utilities > Block
- 2 Select **On** and press **ACT**. The Block feature is now turned **On**. The HOME screen appears with an open circle to show that you have a special feature turned on.

To turn the Block feature off:

- 1 Open the BLOCK OPTION screen.
MAIN MENU > Utilities > Block
- 2 Select **Off** and press **ACT**. The Block feature is now turned **Off**. The HOME screen appears without the open circle.

The Lock Keypad feature

Lock Keypad prevents accidental keypad presses. When the keypad is locked, you can still press **ESC** to view the STATUS screen, the SENSOR STATUS screen, and the sensor graphs.

A locked keypad is also automatically unlocked when batteries are inserted and when alarms/alerts occur.

To lock the keypad:

- 1 Open the **Lock Keypad** screen.
MAIN MENU > Utilities > Lock Keypad
- 2 Select **Utilities**, then press **ACT**.

- 3 Select **Lock Keypad**, then press **ACT**.
- 4 Press **ACT** to lock the keypad.

To unlock the keypad:

- 1 Hold down **SHIFT** and press the  button. A message shows that the keypad is now unlocked.

Running the Selftest

Selftest is a safety utility that allows you to make sure that your monitor is working properly. Selftest runs separately from the routine tests that run independently while the monitor operates.

During Selftest, your monitor automatically runs tests that cause the monitor to beep, vibrate and show different screens. If the monitor detects a condition such as low battery, Selftest stops and shows a message explaining why it stopped.

Contact the HelpLine or your local representative if any of the tests do not run as described below.

- 1 Open the UTILITIES MENU.
MAIN MENU > Utilities
- 2 Select **Selftest**.
- 3 Press **ACT**. The system automatically runs the Selftest in the following order:
 - a. Screen Test:
The screen appears all black.
 - b. Selftest:
The monitor counts down from 10 (it skips the number 4).
 - c. Tone Test:
You should hear beeps.
 - d. Vibrate Test:
You will feel vibrations.
- 4 After Selftest is finished, TEST COMPLETE shows on the screen. The screen returns to the UTILITIES MENU, then to the HOME screen.

User settings

This feature allows you to save your current settings, restore previously saved settings, and clear your current settings. You can also check the history of your use of this feature, to include: the date and time that you used the feature; and if you saved, restored, or cleared the settings.

The Save Settings feature lets you save a group of system settings. After the settings group is saved, you can restore it at a later time as needed. When you restore a settings group, those settings become the active, current settings. Using the Save Settings and Restore Settings features allows you to easily change your current settings as needed without having to re-enter them.

When you clear your system settings, all your current settings are reset to the default settings (the settings that came with your system). After you clear your current settings, do one of the following before you can begin using your system:

- Restore previously saved settings using the Restore Settings feature
- Manually enter your settings into the monitor

Saving your settings

Use this procedure to save your current settings.

- 1 Open the **UTILITIES MENU**.
MAIN MENU > Utilities
- 2 Select **User Settings**. Hold down **SHIFT** and press **ACT**. The **USER SETTINGS** screen shows with **Save Settings** selected.
- 3 Press **ACT**. If you have saved your settings before, then a message shows the last time you saved your settings and asks if you want to overwrite them.
- 4 Do one of the following:
 - Press **ACT** to save your current settings. The **SETTINGS SAVED** message confirms that your current monitor settings have been saved.
 - Press **ESC** if you do not want to save your current settings.
- 5 To return to the **HOME** screen, press **ESC** until it appears.

Restoring your settings

This procedure shows you how to restore the settings you saved most recently (you must have saved your monitor settings in order to restore them).

- 1 Open the UTILITIES MENU.
MAIN MENU > Utilities
- 2 Select **User Settings**.
- 3 Hold down **SHIFT** and press **ACT**. The USER SETTINGS screen appears.
- 4 Select **Restore Settings**, then press **ACT**. A message shows your options.
- 5 Do one of the following:
 - To restore the settings you saved on the date shown, press **ACT**. The SETTINGS RESTORED message confirms that your settings have been restored.
 - If you do not want to restore the settings, press **ESC**.
- 6 Press **ESC** until the HOME screen appears.
- 7 Check your monitor settings to make sure that your settings were restored correctly. See the *Viewing your current settings* section in the *Using your system* chapter for the steps to use to check your settings.

Clearing your settings

Use this utility only if you want to clear your current monitor settings and restore the system default settings.

CAUTION: Do not clear your monitor settings unless directed by your healthcare professional or by a Medtronic representative. If you clear your monitor settings, then you must re-enter them manually or restore a previously saved version of your settings.

To clear your monitor settings:

- 1 Open the UTILITIES MENU.
MAIN MENU > Utilities
- 2 Select **User Settings**.
- 3 Hold down **SHIFT** and press **ACT**. The USER SETTINGS screen appears.
- 4 Select **Clear Settings**, then press **ACT**. The CONFIRM message asks you to confirm that you want to clear your current settings.
- 5 Do one of the following:
 - If you do not want to clear your current settings, press **ESC**.
 - To clear your current settings, select **Yes**, then press **ACT**.
- 6 At this time, the RESET alarm automatically reminds you to reset your system.

- 7 Press  to view the RESET alarm message.
- 8 Press **ESC**, then press **ACT** to clear the RESET alarm. Several screens appear (this happens automatically when you clear a RESET alarm), then the TIME/DATE SETUP screen appears.
- 9 Reset the time and date. You cannot use the system again until you reset the time and date.
- 10 After you reset the time and date, you must enter settings before you can use the system. You can manually re-enter settings, or you can restore the last settings you saved.
- 11 Do one of the following:
 - If you are using a new sensor, then turn on and start the sensor.
 - If you are using the same sensor, then do a **Find Lost Sensor** to re-establish communication.

Viewing your User Settings history

This feature allows you to view your most recent User Settings operations. The system can store up to 20 User Settings operations for review. There are three types:

- **Save** — Saves your current settings.
- **Restore** — Restores previously saved settings.
- **Clear** — Clears your current settings.

To view the User Settings history:

- 1 Open the UTILITIES MENU.
MAIN MENU > Utilities
- 2 Select **User Settings**.
- 3 Hold down **SHIFT** and press **ACT**. The USER SETTINGS screen appears.
- 4 Select **History**, then press **ACT**. The SETTINGS HISTORY screen shows a list of your most recent User Settings operations, including the type of User Settings operation you performed, and the date and time you did it.
- 5 To view more history information, press .
- 6 To return to the HOME screen, press **ESC** until it appears.

Therapy software

Medtronic provides software to optimize your monitor use. Contact our 24 Hour HelpLine or your local representative for more information.

CareLink Personal™

Your Guardian REAL-Time system provides real-time glucose values so that you can track glucose concentration patterns and possibly identify episodes of low and high blood glucose. It also stores glucose data so that it can be analyzed to track patterns and/or downloaded to CareLink Personal for analysis of historical glucose values.

CareLink Personal is a web-based system designed to help you manage your diabetes. It has many key features:

- It copies (uploads) data from your devices: monitor and supported blood glucose meters.
- Device data is stored on a networked database.
- There is an online logbook where you can record self-reported information, such as how many carbohydrates you consumed
- Uploaded data and other information stored on the system can be viewed through several different types of treatment reports.
- Access to data and personal information on this system is secure.

Reports

Reports show you data gathered from your device(s), and from what you directly entered into the system. Examples of the report data you can get include:

- Monitor operation
- Insulin delivery rates, types and times.
- Blood glucose reading levels and times.
- Sensor glucose reading levels and times.

Reports are displayed in PDF format. They can be viewed online, saved, or printed.

You can share your reports with your healthcare professional. These reports can help you and your healthcare professional to discover trends and other information. This can lead to improved therapy management for greater control.

Logbook

The logbook allows you to enter the following self-reported information:

- Carbohydrates consumed
- Exercise activity

The data from your logbook is used in the reports you create. You can see how, for example, your carbohydrate intake compares to blood glucose levels for the same day or time. You can also use it as an online diary of your diabetes management.

Troubleshooting and alarms/alerts

In this chapter:

- About alarms and alerts
- Why alarms and alerts are important
- What to do when you get an alarm/alert
- Sensor alerts
- Viewing your sensor alert history
- Understanding your sensor alert history
- System alarms
- Viewing your system alarm history
- Understanding your system alarm history
- Troubleshooting the MiniLink transmitter and battery

About alarms and alerts

To receive the full benefits of your Guardian REAL-Time CGM System, all parts of the system must work properly. The system also needs to be calibrated regularly so that your sensor glucose (SG) measurements are as close as possible to your meter blood glucose (BG) measurements.

The system is designed to tell you whenever something happens that could affect your diabetes management. The system uses alarms and alerts that are separated into the following categories:

- Calibration
- Glucose status
- Part replacement
- Communication

The STATUS screen shows the last alarm/alert that you cleared. If there was a low battery alarm, this condition will also display in the STATUS screen for as long as the battery works or until you change the battery.

What to do when you get an alarm/alert

When your system generates an alarm or alert:

- an alarm icon  or alert icon  shows on the screen;
- an alarm or alert message shows on the screen; then
- the monitor changes to the HOME screen.

Alarm/alert beeps gradually get louder and louder until you clear them. If the vibrate mode is on, all alarms/alerts start as vibrations, then change to beeps. The beeps change to a siren if you do not respond within ten minutes, then the siren and/or vibrate repeats every one minute until the alarm/alert is cleared.

Follow these steps when you get an alarm or alert:

To view the alarm/alert:

- 1 From the HOME screen, review the alarm/alert by pressing any button to see the message.
- 2 **Read the entire alarm or alert message. The message includes instructions on how to fix the alarm or alert condition.** Press  to scroll if the message is bigger than your screen.
- 3 Press **ESC** to clear the alarm/alert, then press **ACT** after you read the alarm/alert instructions. The HOME screen appears.
- 4 **Follow the instructions** in the alarm/alert message to fix the alarm/alert condition.
- 5 **Check your settings** (for example, time and date and Glucose Limits) to make sure that they are correct.

Sensor alerts

You may encounter the alerts listed below while using the system. These alerts include information about how to resolve the condition.

WEAK SIGNAL

Reason: Occurs when the monitor does not receive data from the transmitter for a pre-defined period of time (as set in Missed Data).

Action: Move the monitor closer to the transmitter or move the transmitter and the monitor to a new location on your body.

LOST SENSOR

Reason: The monitor has not received a signal from the transmitter for more than 40 minutes.

Action: Make sure that the sensor is inserted properly and that the transmitter and sensor are connected. **Do not disconnect.** Touch the inserted sensor at the back of the assembly to prevent movement and push the transmitter firmly:

- If you hear a click, wait 20 seconds to see a green light on the transmitter flash for 10 seconds to confirm a good connection. The alert was due to the sensor and transmitter not being connected.
- If you do not hear a click when you check the connection, the alert was due to a transmission problem. Move the monitor closer to the transmitter.
- If you hear a click, but do not see a green light on the transmitter, check to make sure that the transmitter is charged.

To find your sensor, use the **Find Lost Sensor** function:

MAIN MENU > Sensor > Sensor Start > Find Lost Sensor

LOW TRANSMTR

Reason: Occurs when the transmitter battery is close to running out of power. The transmitter will continue sending sensor signals for several hours and may last for several days until the battery is depleted.

Action: Recharge your transmitter.

CHARGE TRANSMTR

Reason: The transmitter battery is depleted.

Action: Recharge your transmitter immediately. If you are not using a MiniLink transmitter, then your transmitter cannot be recharged and must be replaced.

CHANGE SENSOR

Reason: If this alert results from two Cal Errors in a row, the sensor needs to be replaced. If this alert happens without the Cal Errors, then the transmitter may not be working properly.

Action: If the alert resulted from two Cal Errors in a row, replace the sensor. If the alert happened without Cal Errors, use the transmitter's tester to make sure your MiniLink transmitter is working properly. Call the HelpLine or your local representative for further assistance.

SENSOR END

Reason: The sensor has reached the end of its life.

Action: Replace your sensor. The sensor has a maximum life of 72 hours (three days). The 72-hour life span of the sensor begins when the monitor receives the first METER BG NOW alert.

CAL ERROR

Reason: An error occurred when entering a new meter BG measurement to calibrate the system. Some possible causes are:

- An incorrect blood glucose measurement was entered from the meter into the monitor.
- The entered blood glucose measurement was not current.
- Your blood glucose is rising or falling rapidly.
- The sensor needs more time to stabilize after being inserted.
- The sensor is no longer reading the sensor glucose correctly.

Action: If you receive a Cal Error, wait until your BG is stable to enter a new meter BG measurement for calibration. If you receive a Cal Error on your second calibration, a CHANGE SENSOR alert will occur. Call our 24 Hour HelpLine or your local representative if you have questions.

METER BG NOW

Reason: A meter BG measurement is needed immediately to calibrate the system so that you can continue receiving sensor glucose measurements.

Action: Take and enter a meter blood glucose measurement.

METER BG BY

Reason: A meter BG measurement is required for calibration purposes by the time that is shown on your monitor. The METER BG BY alert is also known as the Cal Reminder alert.

Action: Enter a blood glucose measurement.

LOW SG

Reason: The sensor glucose number is equal to or lower than your low glucose limit. If you do not set a Low Glucose Limit and do not set the glucose alerts ON, then you will not get a Low SG alert.

HIGH SG

Reason: The glucose number is higher than or equal to your high glucose limit. If you do not set a High Glucose Limit and do not set the glucose alerts ON, then you will not get a High SG alert.

LOW PREDICTED

Reason: The sensor glucose measurements may **reach or go below** your Low Glucose Limit in the length of time you selected for the Low Predictive alert.

HIGH PREDICTED

Reason: The sensor glucose measurements may **reach or go above** your High Glucose Limit in the length of time you selected for the High Predictive alert.

FALL RATE

Reason: The sensor glucose measurements are falling at a rate that is equal to or faster than the Set Fall Rate Limit you selected for the alert.

RISE RATE

Reason: The sensor glucose measurements are rising at a rate that is equal to or faster than the Set Rise Rate Limit you selected for the alert.

SENSOR ERROR

Reason: The sensor signal is either too high or too low.

Action: Clear the alert. Do not change the sensor yet. If the alert repeats, make sure that the sensor is inserted properly, that the sensor and transmitter are connected properly, and that there is no moisture at the connections. If the alert continues, check your MiniLink transmitter with the tester. Call the HelpLine or your local representative if this error repeats.

LOW BATTERY

Reason: The battery is running low on energy and may fail at any time. The backlight does not work when the battery is low. If the alert type is set to Vibrate, then the monitor will change to the audio alert Beep Medium.

Action: If you get this alert, then replace the battery before you go to sleep. Clear (ESC, ACT) this alert, then replace the battery.

Viewing your sensor alert history

The SENSOR ALARM HISTORY screen shows you all of the sensor alerts that have recently occurred (up to a maximum of 36). Each line shows the alert date and time, and the alert type.

To view the sensor alert history:

- 1 Open the SENSOR MENU.

MAIN MENU > Sensor

The SENSOR MENU appears.

- 2 Select **Sensor Alarm History**, then press **ACT**. The SENSOR ALARM HISTORY screen shows recent sensor alerts. The most recent (last) sensor alert is selected and appears at the top of the list. From left to right, the entry shows the alarm/alert: date; time, type (Hi or Lo), and SG value.

The list below shows what each sensor alert displayed on your monitor means.

Sensor Alarm History display = this alert

- CalErr = CAL ERROR
- ChgSen = CHANGE SENSOR
- ChgXtr = CHARGE TRANSMITTER
- DwnSlp = FALL RATE
- Hi = HIGH SG
- Lo = LOW SG
- LoBat = LOW BATTERY
- LowXtr = LOW TRANSMITTER
- mBGby = METER BG BY hh:mm
- mBGnow = METER BG NOW
- PrdHi = HIGH PREDICTED
- PrdLow = LOW PREDICTED
- SenEnd = SENSOR END
- SenErr = SENSOR ERROR
- SenLos = LOST SENSOR
- UpSlp = RISE RATE
- WeakSn = WEAK SIGNAL

System alarms

Alarms put the monitor in Attention mode. You may encounter the alarms listed below while using the system. These alarms include information about how to resolve the alarm condition.

A (Alarm) (nn)

Shows a letter A followed by two numbers that show the alarm number. A-alarms cause the monitor to reset. If this alarm repeats often, then call the HelpLine or your local representative for assistance.

BATT OUT LIMIT

Occurs if the battery is removed from the monitor and more than ten minutes passes before a battery is inserted. This alarm may happen when you remove and reinsert a monitor battery.

Go to the UTILITIES menu and reset the date and time. Then follow the **Find Lost Sensor** procedure to re-establish communication with the sensor.

BUTTON ERROR

Occurs if a button has been continually pressed for more than three minutes.

CHECK SETTINGS

When this alarm is active, check and/or reprogram your settings as necessary (including the time and date).

E (Error) (nn)

Shows a letter E followed by two numbers that show the alarm number. E-alarms reset the monitor and clear all of your monitor settings. You must reset your monitor settings after an E-alarm occurs. You must also turn on the sensor, and do the **Find Lost Sensor** procedure to re-establish communication. If you need help, note the error number, then call the HelpLine or your local representative for assistance.

FAILED BATT TEST

The battery does not have enough energy to make the monitor work and the battery cannot be used. Install another battery, and always make sure that you install a new battery into the monitor.

OFF NO POWER

The battery has no more energy. Replace the battery immediately. Follow the directions on the screen. Check to make sure that the time is correct on the screen. Reset the time if necessary.

RESET

Occurs when your settings are cleared for one of the following reasons:

- Your settings were cleared using the Clear Settings feature, and the settings have not been entered again.
- A download attempt from the PC is incomplete (applicable only to the optional Medtronic CareLink Therapy Management software, not available in all locations. Check with your healthcare professional).

WEAK BATTERY

May occur if you install a battery that is less than full strength. The monitor will operate normally, but the battery life will be shorter than normal. Always make sure to install a NEW battery in the monitor.

Viewing your system alarm history

The system allows you to view information on previous system alarms/alerts that you have received, including the name of the alarm, and the date and time it occurred. This is done using the ALARM HISTORY screen. You can view up to 36 of your most recent alarms.

- 1 Open the ALARM HISTORY screen.
MAIN MENU > Utilities > Alarm > Alarm History
- 2 Press  to view previous alarms/alerts.
- 3 To view further information about an alarm/alerts, select the alarm/alert and press **ACT**. The details for that alarm or alert show on the screen.

- 4 Press **ESC** to return to the ALARM HISTORY screen.
- 5 You can continue viewing alarms/alerts by repeating step 2 through step 4 above, or return to the HOME screen by pressing **ESC** until it appears.

The list below shows what each sensor alert displayed on your monitor means.


Alarm History display = this alarm

- Axx = Any A-Alarm
- B-Err = BUTTON ERROR
- B-Out = BATTERY OUT LIMIT
- BadBat = FAILED BATT TEST
- ChkSet = CHECK SETTINGS
- Exx = Any E-Alarm
- NoPwr = OFF NO POWER
- Reset = RESET
- WeakBt = WEAK BATTERY

Troubleshooting the MiniLink transmitter and charger


Checking battery status of the charger and transmitter

Question: Why did the flashing green charger light turn off and a longer flashing red charger light turn on during charging?

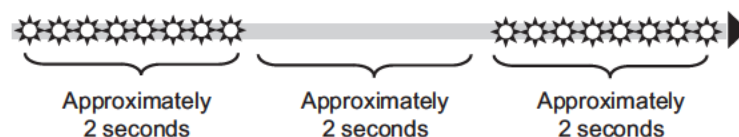


Answer: The transmitter battery is very low. Leave the transmitter on the charger for 8 hours to completely recharge. If the red light still flashes after 8 hours, call the HelpLine or your local representative. It may be time to replace your MiniLink transmitter.

Question: Why do I see quick flashing red lights on the charger?

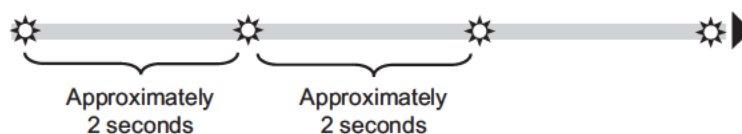

Answer: Your charger battery is low. Replace it with a new AAA battery.

Question: Why do I see a mix of quick and long flashing red lights on the charger?

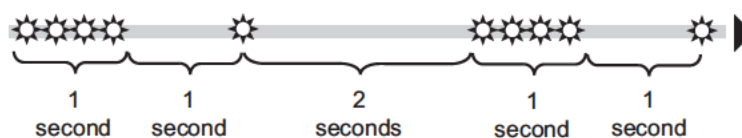

Answer: Your charger battery and MiniLink transmitter batteries are very low. Replace the charger's AAA battery. If you now get the pattern for very low MiniLink transmitter battery, leave the transmitter on the charger for 8 hours to recharge. If the red light is still flashing after 8 hours, call the HelpLine or your local representative. It may be time to replace your MiniLink transmitter.

Question: I had my transmitter on the charger for a day. Will this damage my transmitter?

Answer: It will not damage the transmitter. You cannot overcharge it.

Question: What should I do if the transmitter green light did not flash after removing it from the charger?

Answer: Reconnect the transmitter to the charger for at least one minute, then remove it and watch for the transmitter green light to flash and then turn off.

Question: What should I do if the transmitter green light does not flash when connected to the sensor?

Answer: Make sure the sensor is inserted in your body. If **it is not inserted**, the transmitter will not flash green or send signals to a monitoring device.

If the **sensor is inserted** in your body, disconnect the transmitter from the sensor, wait several seconds and then reconnect it. If the green light still does not flash, charge the transmitter.

Question: Why did I not see the transmitter green light flash after connecting it to the tester?


Answer: Check the connection. If you still do not see a green light flash, fully recharge the transmitter battery. Test the transmitter with the tester. If you still do not see a green light flash, call the HelpLine or your local representative. It may be time to replace your transmitter.

Troubleshooting the transmitter

The blue tester acts like a sensor. If you get sensor alerts, use the tester to make sure that the transmitter is working. If you need help using the tester, contact the HelpLine or your local representative.

Connecting the tester


- 1 Hold the transmitter and the tester as shown. Line up the flat side of the tester with the flat side of the transmitter.

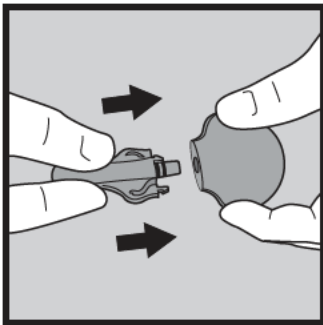

- 2 Push the tester into the transmitter until the flexible side arms of the tester click into the notches on both sides of the transmitter. Within 20 seconds, the green light on the transmitter will flash for approximately 10 seconds when properly connected. You should also see the sensor icon  on your monitor. This means that the transmitter and the sensor are communicating (your monitor Sensor feature must be turned On).
- 3 Go to the RECONNECT OLD SENSOR screen to start your sensor.
Main > Sensor > Sensor Start > Reconnect Old Sensor
- 4 Within the next seven to eight minutes, go to the Sensor Status screen and check the ISIG values:
 - Sensor ISIG value between 24.00 - 29.00 nA - this range of ISIG values means that the transmitter electronics are working properly. Because the transmitter is sending signals correctly, the sensor must have caused the alarm. Remove and discard the sensor. Insert a new sensor in a new site.
 - Sensor ISIG value that is less than 24.00 nA or more than 29.00 nA - call the HelpLine or your local representative. It may be time to replace your transmitter.

Disconnecting the tester

- 1 Hold the transmitter body as shown and pinch the two side arms of the tester.

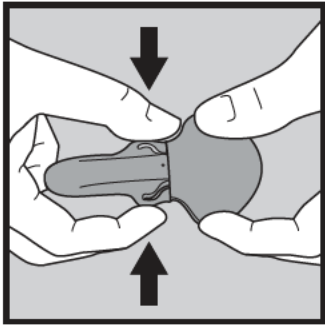

- 2 Pull the transmitter gently from the tester.
- 3 To save transmitter battery energy, do NOT leave the tester connected once you know that the transmitter works.

Troubleshooting the monitor

My monitor will not display my BG measurement from my meter

- 1 Make sure you are using the correct meter (a blood glucose meter supported by MWT1 technology). Your monitor will communicate with this meter only.
- 2 Make sure your meter's RF (radio frequency) feature is on and working correctly.
- 3 Make sure the meter option in your monitor is set to On and you have set the meter ID number correctly on the monitor.
- 4 Make sure your monitor does not have a low battery alert condition.
- 5 Make sure the meter is within 4 feet (1.2 meters) of the monitor without anything in between such as another person, a wall, etc.
- 6 Make sure there is no RF (radio frequency) interference from other electronic devices that could prevent communication. These devices can include some cell phones, cordless phones, televisions, computers, radios, Paradigm pumps, meters and pump remote controls. To restore communication, simply move away from these other devices, or turn them off.

- 7 Your monitor will not show another reading. Make sure the monitor is idle and the HOME screen is blank.
- 8 If your monitor still does not receive your blood glucose reading from the meter, use the up/down buttons to manually enter your blood glucose (in the ENTER BG screen).

Reconnect old sensor

- 1 Go to the RECONNECT OLD SENSOR screen.
Main > Sensor > Sensor Start > Reconnect Old Sensor
- 2 Press **ACT** in the RECONNECT OLD SENSOR screen.
- 3 The SENSOR READY 2 HRS screen appears. Press any button to continue. An alert will sound in two hours when a meter BG is required for calibration.

Find lost sensor

Use this feature if you receive a LOST SENSOR alert, or if communication between the transmitter and the monitor is lost.

- 1 From the HOME screen, open the SENSOR START MENU.
MAIN MENU > Sensor > Sensor Start
The SENSOR START MENU appears with **New Sensor** selected.
- 2 Select **Find Lost Sensor**.
Press **ACT**. A message tells you that your sensor will be ready in 15 minutes.

What happens if I leave the monitor battery out for more than ten minutes?

If you leave the battery out for more than 10 minutes, you may receive a BATT OUT LIMIT alarm message when you install the new battery. If that happens, then do the following:

- 1 Set your monitor time to the correct time, date, and year.
- 2 Make sure that all your settings are correct. If necessary, re-apply your most recently saved settings to the monitor by using the Restore Settings option under User Settings in the UTILITIES MENU (you can only use this option if you have previously saved your monitor settings).
- 3 Check the ALARM HISTORY screen and the STATUS screen for any alarms and/or alerts that may still need attention.
- 4 Perform the **Find Lost Sensor** procedure.

Why doesn't my monitor battery last very long?

- A short battery life does not always mean something is wrong with your monitor. Battery life in your monitor varies, based on the conditions below:
 - The brand of battery you use (Energizer is recommended).
 - The way the battery was stored and/or handled before use (avoid high or low temperatures).
 - How you use your monitor. For example: how often the buttons are pushed, the number of alarms/alerts that occur, and if you keep your graphs on continuous display.
 - Using certain system features (the backlight, vibrate, and meter options) shorten battery life.

What is a CHECK SETTINGS alarm?

This alarm reminds you to make sure that all your settings are correct. A CHECK SETTINGS alarm occurs after:

- The user settings were cleared automatically by the system (set back to their defaults) because there was an E-error alarm.
- You cleared your settings manually using the Clear Settings feature.

My screen looks distorted

The screen may look distorted or have a rainbow-like appearance if you wear polarized sunglasses, are in bright sunlight, or when it is very hot or very cold. This is normal for this type of screen on any device. If your screen looks distorted:

- Take off your sunglasses.
- Move into the shade.
- Make sure that your monitor is not in direct heat (next to a heater) or cold (worn outside your clothing on a very cold day).

I dropped my monitor

- 1 Check the LCD, keypad and monitor case for cracks or damage.
- 2 Review the STATUS screen and your monitor settings to make sure they are still correct.
- 3 Run the Selftest to make sure that the system still works properly.
- 4 Call the HelpLine or your local representative for assistance if necessary.

I submerged my monitor in water

Your monitor is designed to resist brief, accidental contact with water. Do not submerge the monitor in water when bathing, swimming, or other doing water activities.

- 1 Gently dry the outside of the case with a paper or cloth towel. Do not use hot air to dry your monitor. This may damage your monitor's internal electronics.
- 2 Check the battery compartment and the battery. If one or both are wet, allow them to dry completely before using the monitor.
- 3 Run the Selftest to make sure that the system still works properly.

I cannot get to the User Settings screen

If you press **ACT** when you have User Settings highlighted, the following message appears: THIS FEATURE IS NOT NORMALLY ACCESSIBLE. TO ACCESS, CONSULT YOUR USER GUIDE.

To get to the User Settings screen, do the following:

- 1 Go to the UTILITIES MENU.
Main Menu > Utilities
- 2 Select **User Settings**. Hold down **SHIFT** and press **ACT**.
- 3 See the User Settings section of the Utilities chapter for information about the menu options.

System maintenance

In this chapter:

- Replacing the battery
- Storing the system
- Cleaning the system

Replacing the battery

You must replace the monitor battery with a new battery as part of your normal system maintenance. See the Installing the battery section of the Basics chapter for the steps you must follow to replace the monitor battery. These steps include what to do if you receive an alarm after you insert the new battery.

Storing the system

This section shows you how to store the monitor, transmitter, sensor, and Sen-serter.

Storing the monitor

If you must remove and store your monitor for short periods of time, do not remove the battery. For longer-term storage, removing the battery is recommended.

Storing the MiniLink transmitter

If you do not plan to use your MiniLink transmitter for 30 days or more, connect it to the charger for storage. Do not store your MiniLink transmitter in the tester.

Storing the Sen-serter

Store the Sen-serter in the released position to make sure that it lasts as long as possible and works properly.

Cleaning the system

This procedure shows you how to clean the different parts of your system. To avoid accidentally damaging your system during cleaning, be sure to read the Caution section below.

CAUTION:

- Only use a damp cloth and mild liquid soap to clean your system.
 - Never use solvents such as lighter fluid, nail polish remover, or paint thinner to clean your system.
 - Never use lubricants of any kind to clean your system.
-

Before you begin

Make sure that you have the following items before you clean your monitor:

- Clean cloths
- Clean water
- Mixture of clean water and mild liquid soap
- 70 percent alcohol wipes
- Antibacterial hand sanitizer

Cleaning the monitor

This procedure shows you how to clean your monitor.

To clean your monitor:

- 1 Wash your hands thoroughly.
- 2 Dampen a clean cloth with mild liquid soap and warm water, and wipe the outside of the monitor.
- 3 Dampen a clean cloth with clean water.
- 4 Rinse your monitor by wiping it with the dampened part of the cloth.
- 5 Dry your monitor by wiping it with the dry part of the cloth.
- 6 Disinfect your monitor by wiping it with 70 percent alcohol wipes.
- 7 Place the monitor on a clean dry cloth and allow to air-dry for two to three minutes.

CAUTION:

- Never place your monitor directly in water.
 - Keep the battery compartment dry and away from moisture.
 - Never try to clean the inside of your monitor.
-

Cleaning the Sen-serter

This procedure shows you how to clean your Sen-serter.

- 1 Wash your hands thoroughly.
- 2 Dampen a clean cloth with mild liquid soap and warm water, and wipe the Sen-serter.
- 3 Rinse with warm tap water.
- 4 Using an antibacterial hand sanitizer (available at a local drugstore), wipe down the Sen-serter.
- 5 Place the Sen-serter on a clean dry cloth and allow to air-dry for two to three minutes.
- 6 Store the Sen-serter in the released position to make sure that it lasts as long as possible and works properly.

Cleaning the transmitter

CAUTION: The charger and the tester are NOT watertight. Do NOT immerse in water. Do NOT throw away the MiniLink transmitter in a medical waste container or otherwise allow it to be burned. The transmitter contains a battery that may explode when burned.

For single-patient use

- 1 Wash your hands thoroughly.
- 2 Inspect the connector for evidence of blood. If there is any blood inside the connector, the MiniLink transmitter must be discarded and disposed of according to your local laws for non-burning battery disposal.
- 3 Attach the tester to the transmitter to make sure that neither water nor antibacterial hand sanitizer gets inside the connector. Repeated exposure to the hand sanitizer could damage the connectors and affect the transmitter's performance as a result. If you get hand sanitizer inside the connector, allow it to air dry.
- 4 Dampen a clean cloth with mild liquid soap and warm water. Wipe the outside of the transmitter.
- 5 Rinse the transmitter under warm tap water but do not get water inside the connector. If you get water inside the connector, shake the water out and allow it to air dry.
- 6 Using an antibacterial hand sanitizer (available at a local drugstore) on a clean, dry cloth, wipe the transmitter's surface.
- 7 Disconnect the tester from the transmitter.
- 8 Place the transmitter on a clean, dry cloth and air-dry for two to three minutes.

Cleaning the transmitter

WARNING: Thoroughly clean the MiniLink transmitter after removing it from the sensor and before attaching it to the charger. If blood comes in contact with any surface area of the charger, discard the charger immediately. Do not discard the charger in a medical waste container, as it contains a battery that may explode if incinerated. Dispose of the charger according to local regulations for non-incineration battery disposal.

For multiple-patient use

- 1 Wash your hands thoroughly.
- 2 Inspect the connector for evidence of blood. If there is any blood inside the connector, the MiniLink transmitter must be discarded and disposed of according to your local laws for non-burning battery disposal.
- 3 Attach the tester to the transmitter to make sure that neither water nor antibacterial hand sanitizer gets inside the connector. Repeated exposure to the hand sanitizer could damage the connectors and affect the transmitter’s performance as a result. If you get hand sanitizer inside the connector, allow it to air dry.
- 4 Dampen a clean cloth with a mild liquid soap solution. Wipe the outside of the transmitter.
- 5 Rinse the transmitter under warm tap water but do not get water inside the connector. If you get water inside the connector, shake the water out and allow it to air dry.
- 6 Apply three to four drops of a quaternary ammonium compound disinfectant (for example, Cavicide®) on a clean, dry cloth and wipe the transmitter.
- 7 Hold the transmitter and rinse with 70 percent isopropyl alcohol.
- 8 Disconnect the tester from the transmitter.
- 9 Place the transmitter on a clean, dry, non-shedding cloth and air dry.
- 10 After the transmitter is completely dry, place it in a sealed bag labeled with the cleaning date (see label example).

Device:
Date:
Method of Decontamination:
Wash:
Disinfectant:

System specifications

In this chapter:

- Alarm/alert types
- Alarm/alert messages and codes
- Alarm/alert history
- Backlight
- Screen and menu timeout
- Sensor glucose graph timeout
- High and Low glucose limits
- Daily totals
- Default screen
- Default settings
- ENTER METER BG values
- Meter ID entries
- Power supply
- System safety checks
- Monitor size
- Monitor weight
- TIME/DATE SET screen
- STATUS screen
- Guidance and manufacturers declaration

Alarm/alert types

- **Audible tone (beep)** – To make sure that you are aware of the alarm/alert, the sound goes from a beep to a siren until you respond to the alarm/alert.
- **Vibration** – The system vibrates.

Alarm/alert messages and codes

- Alarm/alert information appears on the monitor screen
- Types of alarm/alert information:
 - alarm/alert name
 - alarm/alert message (what to do when you receive the alarm)
 - alarm/alert code (not for all alarms/alerts)

Alarm/alert history

- Maximum number of recent alarms recorded: 36

Backlight

- Type of display: LCD (Liquid Crystal Display)

Screen and menu timeout

- Default: 30 seconds

Sensor glucose graph timeout

- Default: 2 minutes
- Options
 - 2 minutes
 - 4 minutes
 - 6 minutes
 - NONE (no timeout)

High and low Glucose Limits (CSS7100)

- Maximum number of Glucose Limits: 8
- Low Glucose Limit
 - Default 60 mg/dL (3.3 mmol/L)
 - Range 40 to 390 mg/dL (2.2 to 21.6 mmol/L)
 - Increments 1 mg/dL (0.1 mmol/L)
- High Glucose Limit
 - Default 200 mg/dL (11.1 mmol/L)
 - Range 50 to 400 mg/dL (2.8 to 22.2 mmol/L)
 - Increments 1 mg/dL (0.1 mmol/L)

High and low Glucose Limits (CSS7100K)

- Maximum number of Glucose Limits: 8
- Low Glucose Limit
 - Default 90 mg/dL (5.0 mmol/L)
 - Range 90 to 390 mg/dL (5.0 to 21.6 mmol/L)
 - Increments 1 mg/dL (0.1 mmol/L)
- High Glucose Limit
 - Default 280 mg/dL (15.6 mmol/L)
 - Range 100 to 400 mg/dL (5.60 to 22.2 mmol/L)
 - Increments 1 mg/dL (0.1 mmol/L)

Daily totals

- Maximum number of daily totals recorded: 31 days

Default screen

- The HOME screen. The time always shows across the top of the HOME screen.
- The monitor automatically returns to the HOME screen if no buttons are pressed for 30 seconds (for most screens).

Default settings

Feature	Options	Default Setting	Limits, Range of Values
Alarm Clock	• On • Off	Off	Possible values: On Off (You can set up to 8 daily reminders)
Alarm Snooze	(Always on - cannot be turned off)	On 30 minutes	Range: Shortest 0:05 (5 minutes) Longest 1:00 (1 hour)
Alert Silence	• Off • Hi only • Lo only • Hi and Lo • All	Off	Possible values: Off High BG Low BG All
Alert Type	• Beep • Vibrate	Beep Medium	Possible values: Beep Beep Long Beep Medium Beep Short Vibrate

Feature	Options	Default Setting	Limits, Range of Values
AUC Limits	Always on	Low Limit 70 (mg/dL) 3.9 (mmol/L) High Limit 180 (mg/dL) 10.0 (mmol/L)	Low Limit range: ^a Lowest 40 (mg/dL) 2.2 (mmol/L) Highest 400 (mg/dL) 22.2 (mmol/L) High Limit range: ^b Lowest 40 (mg/dL) 2.2 (mmol/L) Highest 400 (mg/dL) 22.2 (mmol/L)
Block	- On - Off	Off	Possible values: On Off
Blood Glucose Units (BG Units)	- mg/dL - mmol/L	mg/dL	Possible values: mg/dL mmol/L
Cal Reminder	- On - Off	On 1 hour	Range: Shortest 0:05 (5 minutes) Longest 6:00 (6 hours)

a. The AUC Low Limit can be equal to, but not higher than the AUC High Limit.

b. The AUC High Limit can be equal to, but not lower than the AUC Low Limit.

Feature	Options	Default Setting	Limits, Range of Values
Glucose Alarms	- On - Off	Off	You set the thresholds (limits) for the Glucose Alarms using the Glucose Limits settings.
Glucose Limits	- Maximum of 8 Glucose Limits. - One pair of Glucose Limits (Glucose Limits 1) is already set up for you, with a default High Glucose Limit and Low Glucose Limit. - Glucose Limits 2 through 8 can be set as needed.		
Glucose Limits 1	CSS7100 Low: Off, On High: Off, On CSS7100K Low - Off, or valid range from 90 to 390 mg/dL (5.0 mmol/L to 21.6 mmol/L) High - Off, or valid range from 100 mg/dL to 400 mg/dL (5.6 mmol/L to 22.2 mmol/L)	CSS7100 Low: On 60 mg/dL (3.3 mmol/L) High: On 200 mg/dL (11.1 mmol/L) CSS7100K Low: On 90 mg/dL (5.0 mmol/L) High: On 280 mg/dL (15.5 mmol/L)	Low Glucose Limit range:^a - CSS7100: 40 - 390 mg/dL (2.2 - 21.6 mmol/L) - CSS7100K: 90 - 390 mg/dL (5.0 - 21.6 mmol/L) High Glucose Limit range:^b - CSS7100: 50 - 400 mg/dL (2.8 - 22.2 mmol/L) - CSS7100K: 100 - 400 mg/dL (5.6 - 22.2 mmol/L)
a. Your Low Glucose Limit must be at least 10 mg/dL (0.6 mmol/L) below your High Glucose Limit. For example, if your High Glucose limit is 200 mg/dL (11.1 mmol/L), then your Low Glucose limit cannot be higher than 190 mg/dL (10.5 mmol/L). b. Your High Glucose Limit must be at least 10 mg/dL (0.6 mmol/L) above your Low Glucose Limit.			

Feature	Options	Default Setting	Limits, Range of Values
Glucose Limits 2 through Glucose Limits 8	CSS7100 Low: Off, On High: Off, On CSS7100K Low - Off, or valid range from 90 to 390 mg/dL (5.0 mmol/L to 21.6 mmol/L) High - Off, or valid range from 100 mg/dL to 400 mg/dL (5.6 mmol/L to 22.2 mmol/L)	CSS7100 Low Off High Off CSS7100K Low Off High Off	See the Glucose Limits 1 range.
High Snooze	- On - Off (Set to On when High Glucose, High Predictive, or Rising Rate of Change Alarm is on)	On 1 hour	Range: Shortest 0:05 (5 minutes) Longest 3:00 (3 hours)
a. Your Low Glucose Limit must be at least 10 mg/dL (0.6 mmol/L) below your High Glucose Limit. For example, if your High Glucose limit is 200 mg/dL (11.1 mmol/L), then your Low Glucose limit cannot be higher than 190 mg/dL (10.5 mmol/L). b. Your High Glucose Limit must be at least 10 mg/dL (0.6 mmol/L) above your Low Glucose Limit.			

Feature	Options	Default Setting	Limits, Range of Values
Language	- Cesky - Dansk - Deutsch - English - Español - Français - Italiano - Magyar - Nederlands - Norsk - Polski - Portugues - Slovenscina - Suomi - Svenska - Turkse	English	Possible values: Cesky Dansk Deutsch English Español Français Italiano Magyar Nederlands Norsk Polski Portugues Slovenscina Suomi Svenska Turkse
Lock Keypad	- On - Off	Off	Possible values: On Off
Low Snooze	- On - Off (Set to On when Low Glucose, Low Predictive, or Falling Rate of Change Alarm is on.)	On 20 minutes	Range: Shortest 0:05 (5 minutes) Longest 1:00 (1 hour)
Meter Option	- On - Off	Off	Possible values: On Off

Feature	Options	Default Setting	Limits, Range of Values
Missed Data	(Always on - cannot be turned off)	30 minutes	Range: Shortest 0:05 (5 minutes) Longest 0:40 (40 minutes)
Predictive Alarms	- On - Off	Off Default time sensitivity: Low 15 minutes High 15 minutes	Range: Off Shortest 0:05 (5 minutes) Longest 0:30 (30 minutes)
Rate of Change Alarms	- On - Off	Fall Rate Limit Off Rise Rate Limit Off	Fall Rate Limit range: Lowest 1.1 (mg/dL/min) 0.065 (mmol/L/min) Highest 5.0 (mg/dL/min) 0.275 (mmol/L/min) Rise Rate Limit range: Lowest 1.1 (mg/dL/min) 0.065 (mmol/L/min) Highest 5.0 (mg/dL/min) 0.275 (mmol/L/min)
Sensor	- On - Off	Off	Possible values: On Off

Feature	Options	Default Setting	Limits, Range of Values
Sensor glucose graph timeout	2 minutes 4 minutes 6 minutes - NONE^a	2 minutes	Possible values: Shortest 0:02 (2 minutes) 0:04 (4 minutes) 0:06 (6 minutes) No timeout NONE
Time/Date set-up	12 Hour 24 Hour	12 Hour	Possible values: 12 Hour 24 Hour
a. If you select this option the sensor glucose graphs will show continuously on the monitor screen.			

ENTER BG meter values

- Expiration: 12 minutes
- Supported value range:
 - from 1.1 to 33.3 (mmol/L)
 - from 20 to 600 (mg/dL)
- Appears in the STATUS screen (RF only)

Meter ID entries

- Maximum meter ID entries: 3

Power supply

- Standard 1.5 V AAA alkaline battery
- Energizer brand is recommended

System safety checks

- More than 50 independent checks
- Regular and periodic checking of system functions

Monitor size

- Inches
 - Height approximately 2.0 inches
 - Width approximately 3.2 inches
 - Depth (thickness) approximately 0.77 inches
- Centimeters (cm)
 - Height approximately 5.1 cm
 - Width approximately 8.1 cm
 - Depth (thickness) approximately 2.0 cm

Monitor weight

- Approximately 79 grams (with battery installed)

TIME/DATE SET screen

- 12-hour or 24-hour formats
- You set the time/date (including the year, month and day)

STATUS screen

Screen item	Description
Last Alarm/alert:	The most recent alarm/alert that occurred within the past 24 hours. Alarms that occurred more than 24 hours ago are not displayed.
SENSOR DEMO: On	Shows if the Sensor Demo feature is turned On or Off . This only appears when the Sensor Demo feature is turned On .
Meter: On	Shows if the meter option is turned On or Off . This only appears when the meter option is turned On .
Meter: On, Low Batt	Shows that the meter option is turned On and the monitor battery is low or out of energy. This only appears when the meter option is turned On and the monitor battery is low or out of energy.
Battery:	Shows battery condition (how much energy is left in the battery). This always appears.
BG meter value:	Shows the most recent meter blood glucose value that you entered within the past 24 hours, and when it was entered (date and time). Meter blood glucose values that were entered more than 24 hours ago are not displayed. This only appears when the meter option is turned on.
Date:	Shows the current date. This always appears.
S/N# XXXXXX	The serial number of the monitor. This always appears.
Guardian XXXX	The model number of the monitor. This always appears.
Software ver- sion:	The version number of the system software. This always appears.
Status of system	Always appears (for example, Set Time).

Guidance and manufacturer's declaration

Guidance and Manufacturer's Declaration - Electromagnetic Emissions		
The Guardian REAL-Time system (consisting of the CSS7100/CSS7100K Receiver and MMT-7703 transmitter) is intended for use in the electromagnetic environment specified below. The customer or the user of the Guardian REAL-Time system should make sure that it is used in such an environment.		
Emissions Test	Compliance	Electromagnetic Environment - Guidance
RF emissions Test: CISPR 11	Group 1	The Guardian REAL-Time system uses RF energy only for system communication functions. Therefore, its RF emissions are very low and are not likely to cause any interference in nearby electronic equipment.
RF emissions CISPR 11 Harmonic emissions IEC 61000-3-2 Voltage fluctuations/flicker emissions IEC 61000-3-3	Class B Not applicable Not applicable	The Guardian REAL-Time system is suitable for use in all establishments other than domestic and those directly connected to the public low-voltage power supply network that supplies buildings used for domestic purposes. NOTE: The preceding statement is required by IEC 60601-1-2 for Group 1, Class B devices. However, since the Guardian REAL-Time system is battery powered, its emissions will not be affected by the establishment power supply and there is no evidence of any issues associated with the use of the system in domestic establishments.

Guidance and Manufacturer's Declaration - Electromagnetic Immunity			
The Guardian REAL-Time system (consisting of the CSS7100/CSS7100K Receiver and MMT-7703 transmitter) is intended for use in the electromagnetic environment specified below. The customer or the user of the Guardian REAL-Time system should make sure that it is used in such an environment.			
Immunity Test	IEC 60601 Test Level	Compliance Level	Electromagnetic Environment - Guidance
Electrostatic discharge (ESD) IEC 61000-4-2	±6 kV contact ±8 kV air	±30 kV air (<5% relative humidity)	The Guardian REAL-Time system should not be affected by electrostatic discharge that might occur under normal conditions of use.
Electrical fast transient/burst IEC 61000-4-4	±2 kV for power supply lines	Not applicable	Not applicable
	±1 kV for input/output lines	Not applicable	Not applicable
Surge IEC 61000-4-5	±1 kV line(s) to line(s)	Not applicable	Not applicable
	±2 kV line(s) to earth	Not applicable	Not applicable
Voltage dips, short interruptions and voltage variations on power supply lines IEC 61000-4-11	<5% U_T (>95% dip in U_T) for 0.5 cycle	Not applicable	Not applicable
	40% U_T (60% dip in U_T) for 5 cycles	Not applicable	Not applicable
	70% U_T (30% dip in U_T) for 25 cycle	Not applicable	Not applicable
	<5% U_T (>95% dip in U_T) for 5 seconds	Not applicable	Not applicable

Guidance and Manufacturer's Declaration - Electromagnetic Immunity			
The Guardian REAL-Time system (consisting of the CSS7100/CSS7100K Receiver and MMT-7703 transmitter) is intended for use in the electromagnetic environment specified below. The customer or the user of the Guardian REAL-Time system should make sure that it is used in such an environment.			
Immunity Test	IEC 60601 Test Level	Compliance Level	Electromagnetic Environment - Guidance
Power frequency (50/60 Hz) magnetic field IEC 61000-4-8	3 A/m	3 A/m	Power frequency magnetic fields should be at levels characteristic of a typical location in a typical commercial or hospital environment.
NOTE: U_T is the a.c. mains voltage prior to application of the test level.			

Guidance and Manufacturer's Declaration - Electromagnetic Immunity			
The Guardian REAL-Time is intended for use in the electromagnetic environment specified below. The customer or user of the Guardian REAL-Time system should make sure that it is used in such an environment.			
Immunity Test	IEC 60601 Level	Compliance Level	Electromagnetic Environment Guidance
Conducted RF IEC 61000-4-6	10Vrms 150kHz to 80MHz	Not applicable	Portable and mobile RF communications equipment should be used no closer to any part of the Guardian REAL-Time, including cables, than the recommended separation distance calculated from the equation applicable to the frequency of the transmitter. Recommended separation distance Not applicable
Radiated RF IEC 61000-4-3	10V/m 80MHz to 6.0GHz	10 V/m	$d = 0.35 \sqrt{P}$ 80 MHz to 800 MHz $d = 0.70 \sqrt{P}$ 800 MHz to 6.0 GHz Where P is the maximum output power rating of the transmitter in watts (W) according to the transmitter manufacturer and d is the recommended separation distance in meters (m). Field strengths from fixed RF transmitters, as determined by an electromagnetic site survey^a, should be less than the compliance level in each frequency range^b. Interference may occur in the vicinity of equipment marked with the following symbol: 

Guidance and Manufacturer's Declaration - Electromagnetic Immunity			
The Guardian REAL-Time is intended for use in the electromagnetic environment specified below. The customer or user of the Guardian REAL-Time system should make sure that it is used in such an environment.			
Immunity Test	IEC 60601 Level	Compliance Level	Electromagnetic Environment Guidance
NOTE: At 80 MHz and 800 MHz, the higher frequency range applies. NOTE: These guidelines may not apply in all situations. Electromagnetic propagation is affected by absorption, and reflection from structures, objects and people. ^aField strengths from fixed transmitters, such as base stations for radio (cellular/cordless) telephones and land mobile radios, amateur radio, AM and FM radio broadcasts and TV broadcast cannot be predicted theoretically with accuracy. To access the electromagnetic environment due to fixed RF transmitters, an electromagnetic site survey should be considered. If the measured field strength in the location in which the Guardian REAL-Time is used exceeds the applicable RF compliance level above, the Guardian REAL-Time should be observed to verify normal operation. If abnormal performance is observed, additional measures may be necessary, such as re-orienting or relocating the Guardian REAL-Time. ^bOver the frequency range 150 kHz to 80 MHz, field strengths should be less than [V₁] V/m.			

Recommended separation distances between portable and mobile RF communications equipment and the Guardian REAL-Time system

This section provides information on the recommended separation distance between portable and mobile RF communications equipment and the Guardian REAL-Time. The Guardian REAL-Time is intended for use in an electromagnetic environment in which radiated RF disturbances are controlled. Guardian REAL-Time users can help prevent electromagnetic interference by maintaining a minimum distance between portable and mobile RF communications equipment (transmitters) and the Guardian REAL-Time as recommended below, according to the maximum output power of the communications equipment.

	Separation distance according to the frequency of transmitter (m)	
Rated maximum output power of transmitter (W)	80MHz to 800MHz $d = 0.35 \sqrt{P}$	800MHz to 6.0GHz $d = 0.70 \sqrt{P}$
0.01	0.035	0.07
0.1	0.11	0.11
1	0.35	0.7
10	1.1	2.2
100	3.5	7

For transmitters rated at a maximum output power not listed above, the recommended separation distance d in meters (m) can be estimated using the equation applicable to the frequency of the transmitter, where p is the maximum output power rating of the transmitter in watts (W) according to the transmitter manufacturer.

NOTE: At 80 MHz and 800 MHz, the separation distance for the higher frequency range applies.

NOTE: These guidelines may not apply in all situations. Electromagnetic propagation is affected by absorption and reflection from structures, objects and people.

MiniLink specifications

Biocompatibility	MiniLink: Complies with ISO 10993-1 for long-term body contact
Operating Conditions	MiniLink Temperature: +32° to +122°F (0° to 50°C) MiniLink Relative Humidity: 10% to 95% with no condensation Charger Temperature: +50° to +104°F (10° to 40°C) Charger Relative Humidity: 30% to 75% with no condensation
Storage Conditions	MiniLink Temperature: -4° to +131°F (-20° to +55°C) MiniLink Relative Humidity: 10% to 100% with no condensation Charger Temperature: +14° to +122°F (-10° to +50°C) Charger Relative Humidity: 10% to 95% with no condensation
Battery Life	MiniLink: 14 days of continuous glucose monitoring immediately following a full charge Charger: Completes 40 typical charging operations using a new AAA battery
Transmitter Frequency	868.35 (WW configuration) or 916.5 megahertz (NA configuration)

Symbol table

Do not reuse	
Attention: See Instructions for Use	
Date of manufacture (year - month)	
Manufacturer	
Batch code	LOT
Use by (year - month)	
Catalogue number	REF
Device serial number	SN
Configuration	CONF
Storage temperature range	
Fragile product	
Type BF equipment (Protection from electrical shock)	
Monitor: Conforms to IEC60601-1 sub-clause 44.6 and IEC60529 standard	IPX7
Transmitter: Protected against the effects of continuous immersion in water (8 feet or 2.4 meters immersion for 30 minutes).	IPX8
Recycle	
Radio communication	

One Per Container /Package	(1X)
CE mark by notified body as a medical device	
CE mark as a radio transmitter under the R & TTE 1999/5/EC/directives. Applies to 868.35 megahertz (WW configuration):	
Signifies compliance with Industry Canada EMC and Radiocommunications requirements:	IC
Signifies compliance with Australian EMC and Radiocommunications requirements:	
Magnetic Resonance (MR) unsafe:	

Sensor accuracy

The information in this section is meant to be reviewed by you with your healthcare professional. Please also note: The Guardian RT is not an abbreviation for the Guardian REAL-Time CGM System. The Guardian RT is a continuous glucose monitoring system that was developed before the Guardian REAL-Time CGM System. The Guardian RT was originally named TGMS II, and the name was changed to Guardian RT.

The Guardian REAL-Time CGM System uses a glucose sensor to continuously monitor your glucose levels. The Guardian REAL-Time CGM System uses the same algorithm as the Guardian RT®. The Guardian RT was evaluated in two clinical studies and the following information in this section describes the findings. The remainder of this section discusses the results of this pivotal study, and uses the term Guardian RT.

The Medtronic Guardian RT uses a glucose sensor to continuously monitor your glucose levels. The Guardian RT sensor is calibrated using your home blood glucose meter. Once calibrated, the Guardian RT reports glucose values every five minutes. These values were compared to reference laboratory blood glucose measurements to check the Guardian RT's performance characteristics in two clinical studies.^{1, 2}

Although presentations to characterize performance of the Guardian RT are given below, there is no commonly accepted statistical approach for capturing the performance of continuous glucose monitors such as the Guardian RT. Performance may be best characterized by viewing graphs called time-elapsed plots. In these plots, the values from Guardian RT for one subject over time are overlaid with values at the same time from the glucose reference method. Three representative time-elapsed plots are shown at the end of this section of sensors that exhibited excellent performance, average performance and poor performance.

1. Medtronic Diabetes, A Frequent Sample Accuracy Evaluation of the Medtronic Diabetes Telemetered Glucose Monitoring System II (TGMS II) in Subjects with Type 1 Diabetes Mellitus, August 2004.
2. Medtronic Diabetes, An Accuracy Evaluation of the Medtronic Diabetes Guardian RT Glucose Monitoring System in Pediatric Subjects with Type 1 Diabetes Mellitus, February, 2006.

Performance results in adults

The performance of the Guardian RT in adults was evaluated in a clinical study.³ Guardian RT results were compared to plasma glucose values from a reference method, the YSI 2300 STAT Plus™ glucose analyzer (referred to as YSI). Sixteen subjects with Type I diabetes participated in a single-site in-clinic study. Subjects ranged in age from 18 to 65 years old. Each subject wore two Guardian RT systems simultaneously. One Guardian system was calibrated an average of 3.5 times per day, and the other was calibrated approximately five times per day using the Paradigm Link blood glucose meter. YSI measurements were taken every 30 minutes.

Users and their healthcare professionals should consider that performance in this study may be idealized, and that performance may be worse when the Guardian RT is used in a less controlled home setting. For example:

- The mean Hemoglobin A1c among the 16 participants was 8.2 percent. As hemoglobin A1c levels rise, conditions often occur which are most challenging to test systems measuring glucose in interstitial fluid, i.e., higher glucose levels, more rapid changes in glucose concentrations, and often more hypoglycemic episodes.
- Subjects saw, on average, between four and five fingerstick values per day. This enables subjects to better manage their diabetes when compared to those who perform less fingersticks per day. Agreement between Guardian RT and YSI values is shown to be closer at mid-range glucose levels, as compared to agreement at low or high glucose concentrations.
- Subjects were more limited in their activities than what may exist in home use, and they were provided with all their meals. Sensors were also inserted by clinic staff rather than the subjects themselves. Subjects who are more active, or with poor eating habits, may create more challenging conditions for the Guardian RT.
- Performance of the Guardian RT may vary depending on the glucose meter used and how well the meter is maintained. It is important to carry out quality-control checks on the meter and code the meter according to the manufacturer's instructions to optimize performance of the Guardian RT.

3. Medtronic Diabetes, A Frequent Sample Accuracy Evaluation of the Medtronic Diabetes Telemetered Glucose Monitoring System II (TGMS II) in Subjects with Type 1 Diabetes Mellitus, August 2004.

Accuracy of Guardian RT readings

In this study, YSI measurements (taken every half hour) were paired with the corresponding Guardian RT reading (taken every five minutes). Pairing was done by selecting the Guardian RT value closest in time to the YSI test result. Agreement was analyzed by comparing paired glucose measurements.

Agreement between the matched pair was estimated by evaluating the difference between the Guardian RT reading and the YSI measurement. The difference between them was calculated as a percentage of the YSI (Mean Absolute Percent Difference). The bias was also calculated, and it is defined as the overall difference between the Guardian RT glucose values and the YSI values. The paired glucose measurements are summarized below.

Number of Paired Glucose Measurements	3941
Mean Absolute Percent Difference ($\pm$ SD)	19.7 $\pm$ 18.4%
Bias	-0.8 mmol/L (-15.0 mg/dL)

The accuracy of the Guardian RT was also evaluated by calculating the percentage of Guardian RT readings within 20 percent and within 30 percent of the YSI reading (or within 1.1 mmol/L (20 mg/dL) in the low glucose range). Results are shown below.

Plasma Glucose Range (mg/dl)	Plasma Glucose Range (mmol/l)	Number of Paired Readings	Percent Within 20%	Percent Within 30%
Overall		3941	62%	79%
40-80 ^a	2.2-4.4	356	68%	68%
>80-120	>4.4-6.7	769	60%	77%
>120-240	>6.7-13.3	2362	62%	81%
>240	>13.3	454	61%	82%

a. For the Low glucose range, 2.2-4.4 mmol/L (40-80 mg/dL), the value shown is the percent within 1.1 mmol/L (20 mg/dL).


The Clarke Error Grid was used to assess the clinical relevance of the differences between the Guardian RT readings and the comparative YSI measurements. The Clarke Error Grid divides a correlation plot into 5 zones. See following table.

Results in zones A and B are considered clinically acceptable, while results in zones C, D, and E are potentially dangerous and, therefore, clinically significant errors. The Clarke Error Grid zones are labeled on the correlation plot.

Zone	Description
A	Clinically accurate, would lead to correct treatment decisions
B	Would lead to benign decisions or no treatment
C	Would lead to over-correction of normal glucose levels
D	Would lead to failure to detect and treat high or low glucose levels
E	Would lead to erroneous treatment decisions

The plot below is a correlation plot of Guardian RT readings versus readings from the reference method, the YSI 2300 Glucose Analyzer. It is overlaid with the Clarke Error Grid.

The total number of paired data points is 3941.

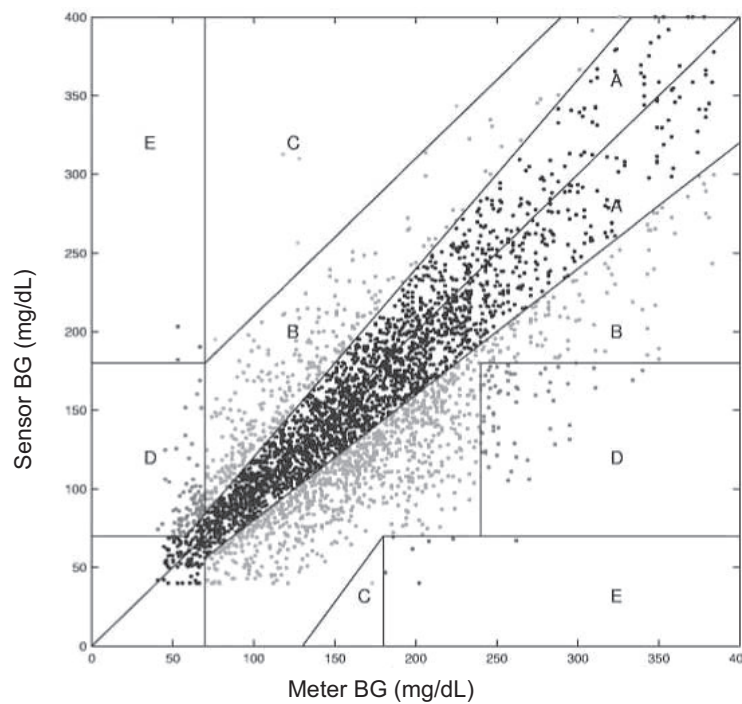


The percent of Guardian RT readings in the previous graph are presented in the following table according to the percentage of points falling within each zone (A-E). Results are further broken down (stratified) according to the range of glucose concentrations.

Glucose Range (mg/dl)	Number and (%) of Data Points Evaluated	A + B	A	B	C	D	E
40-80	356 (9)	271 (76.1)	214 (60.1)	57 (16.0)	2 (0.6)	80 (22.5)	3 (0.8)
81-120	769 (20)	768 (99.9)	463 (60.2)	305 (39.7)	1 (0.1)	N/A ^a	N/A
121-240	2362 (60)	2352 (99.6)	1476 (62.5)	876 (37.1)	4 (0.2)	N/A	6 (0.2)
>240	454 (11)	394 (86.8)	277 (61.0)	117 (25.8)	N/A	59 (13.0)	1 (0.2)
Overall	3941 (100)	3785 (96.0)	2430 (61.7)	1355 (34.4)	7 (0.2)	139 (3.5)	10 (0.2)
a. N/A means that the Clarke Error Grid does not consider the possibility of these zones in that concentration range.							

Precision of Guardian RT readings

This study was also designed to look at the reproducibility of two sensors worn simultaneously at different locations on the body. Precision was estimated by comparing the glucose readings from the two Guardian RT systems. In this study 11,475 paired sensor Guardian RT values were obtained. On average, they were different by 17.2 percent. **The following figure is an example of how data was paired in this study.** In the graph there are two tracings of Guardian RT values. Each tracing comes from a different Guardian RT unit worn by one subject during a one-day period.

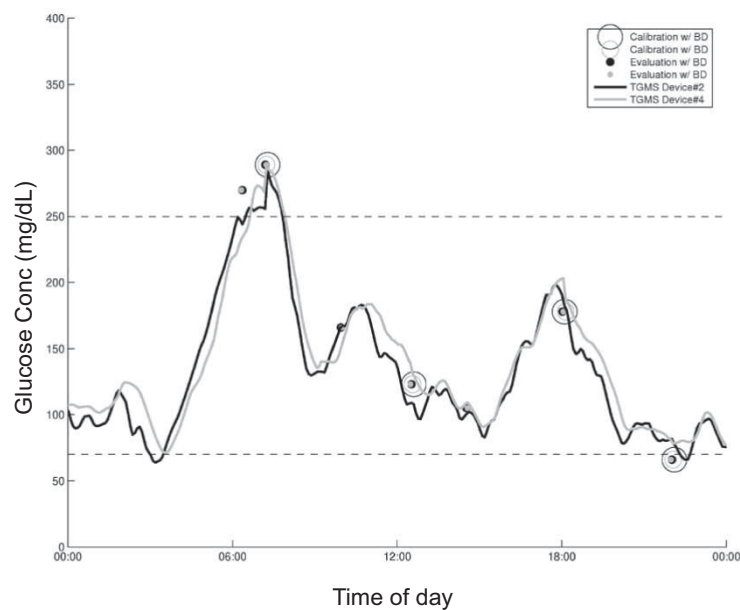

The above chart, Guardian RT device #2 is represented by a solid line, and the Guardian RT device #4 by a light-gray line.

Low and High Alerts in adults

The ability of the Guardian RT to detect high and low glucose levels was measured in the same clinical study.⁶ Since it is important to set the alert levels in a conservative fashion, the Low Glucose Alert should be set at a value slightly higher than the value of blood glucose you want to detect, and the High Glucose Alert should be set at a value slightly lower than the value of blood glucose you want to detect.

NOTE: Please ask your healthcare professional which low and high alert setting is best for you.

The Low Glucose Alert

The Low Glucose Alert was evaluated for its ability to detect glucose levels at 70 mg/dL (3.9 mmol/L), or below, using the YSI 2300 STAT Plus glucose analyzer. As a reference, with the Low Glucose Alert set at 70 mg/dL (3.8 mmol/L), 49% (100/205) of low glucose events were detected by the Guardian RT. Better detection of low blood glucose can be obtained by setting the Low Glucose Alert level higher. For example, setting the Low Glucose Alert at 90 mg/dL (5.0 mmol/L), instead of 70 mg/dL (3.9 mmol/L), increases the ability to detect low blood glucose levels from 49% to 82% (see the table below).

Sometimes the Guardian RT will alert when the blood glucose levels are not low. When the Guardian RT Low Alert was set at 70 mg/dL (3.9 mmol/L) in this study, 43% of the results were considered false alerts (actual blood glucose values are greater than 85 mg/dL (4.7 mmol/L)). This percentage may be exaggerated because blood glucose may be dropping when the Guardian RT alerts. The table below shows the percent of Low Glucose readings correctly identified by the Guardian RT for specific settings.

Guardian RT Low Alert Setting (mg/dL)	Guardian RT Low Alert Setting (mmol/L)	True Alert Rate ^a	False Alert Rate ^b
70	3.9	49%	60%
80	4.4	68%	64%

a. True Alert Rates are the % of times when the glucose level was at or below the alert setting and the alert sounded.

b. False Alerts Rates are the % of times when the Guardian RT Sensor alarmed but the blood glucose level was greater than the alert setting.

6. Medtronic Diabetes, A Frequent Sample Accuracy Evaluation of the Medtronic Diabetes Telemetered Glucose Monitoring System II (TGMS II) in Subjects with Type 1 Diabetes Mellitus, August 2004.

Guardian RT Low Alert Setting (mg/dL)	Guardian RT Low Alert Setting (mmol/L)	True Alert Rate ^a	False Alert Rate ^b
90	5.0	82%	75%
100	5.6	90%	79%

a. True Alert Rates are the % of times when the glucose level was at or below the alert setting and the alert sounded.

b. False Alerts Rates are the % of times when the Guardian RT Sensor alarmed but the blood glucose level was greater than the alert setting.

Increasing the Low Alert settings will improve the ability to detect low blood glucose events, but it will also increase the frequency of Guardian RT false alerts for blood glucose levels not below the target value. You should consider this trade-off between the improved ability to detect true low blood glucose versus the increased number of false alerts when setting the low alert threshold.

The High Glucose Alert

The High Glucose Alert was evaluated for its ability to detect glucose levels at 250 mg/dL (13.8 mmol/L), or above, using the YSI analyzer. As a reference, with the High Glucose Alert set at 250 mg/dL (13.8 mmol/L), 53% (195/365) of high glucose events were detected by the Guardian RT. Better detection of high blood glucose can be obtained by setting the High Glucose Alert level lower. For example, setting the High Glucose Alert at 190 mg/dL (10.6 mmol/L), instead of 250 mg/dL (13.8 mmol/L), increases the ability to detect high blood glucose levels from 53% to 85% (see the table below).

Sometimes the Guardian RT will alert when the blood glucose levels are not high. When the Guardian RT High Alert was set at 250 mg/dL (13.8 mmol/L) in this study, 7.2% of the results were considered false alerts (actual blood glucose values are less than 225 mg/dL (12.5 mmol/L)). This percentage may be exaggerated because blood glucose may be rising when the Guardian RT alerts. The table below shows the percent of High Glucose readings correctly identified by the Guardian RT for specific settings.

Guardian RT High Alert Setting (mg/dL)	Guardian RT High Alert Setting (mmol/L)	True Alert Rate*	False Alert Rate**
190	10.6	85%	64%
200	11.1	81%	58%

Guardian RT High Alert Setting (mg/dL)	Guardian RT High Alert Setting (mmol/L)	True Alert Rate*	False Alert Rate**
225	12.5	67%	40%
250	13.8	53%	25%

* True Alert Rates are the % of times when the glucose level was at or above the alert setting and the alert sounded.

** False Alerts Rates are the % of times when the Guardian RT Sensor alarmed but the blood glucose level was lower than the alert setting.

Decreasing the High Alert settings will improve the ability to detect high blood glucose events, but it will also increase the frequency of Guardian RT false alerts for blood glucose levels not above the target value. You should consider this trade-off between the improved ability to detect true high blood glucose versus the increased number of false alerts when setting the high alert threshold.

Guardian RT Sensor Performance and Calibration Stability As a Function of Time

The Guardian RT sensor may be worn for up to 3 days (72 hours) and must be calibrated at least twice a day. Two sets of data, approximately equal in number, were collected during the clinical trial. One data set was generated when the frequency of calibrations averaged 3.5 per day (Data Set A), and the other averaged 5 times a day (Data Set B). During the study, a total of 38 sensors were evaluated in 16 individuals.

As per the stratified Clarke Error Grid analysis above, agreement between Guardian RT values and YSI values tends to be less at low and high glucose concentrations when compared to other concentration ranges.

Guardian RT performance in the hypoglycemic range, as a function of sensor insertion time, is characterized below. Results from the two different data sets are presented. The two populations were separated according to the number of calibrations per day. This table represents the percentage of Data Points in the 2.2 - 4.4 mmol/L (40-80 mg/dL) range that fell within 1.1 mmol/L (20 mg/dL). Data is presented in 12-hour increments.

Data Set	0-12 hrs	12-24 hrs	24-36 hrs	36-48 hrs	48-60 hrs	60-72 hrs
A	78%	81%	73%	65%	56%	41%
B	67%	70%	93%	60%	75%	38%

An analysis of the mean percentage of Absolute Relative Error (ARE percentage) and standard deviations, across 12-hour increments of wear periods, appears in the table below. Both data sets are pooled together in this data.

Hours From Insertion	Mean ARE (%)	Std. Dev.
0-12 hrs	24.84	20.04
12-24 hrs	19.66	16.17
24-36 hrs	16.43	15.62
36-48 hrs	18.23	19.27
48-60 hrs	16.59	14.25
>60 hrs	22.95	23.51

The median sensor life from Data Sets A and B were 57.5 hours and 72.9 hours, respectively. Twenty-one of the sensors operated for 72 hours, while the others were removed for a variety of reasons, most often because of calibration errors.

The percentage of Guardian RT readings within 20 percent and 30 percent of YSI readings was analyzed according to time after sensor insertion and according to the glucose-concentration range (as determined by the YSI analyzer). See the following table.

Glucose Range (mg/dl)	Percentage of Guardian RT values within 20% of YSI laboratory readings		Percentage of Guardian RT values within 30% of YSI laboratory readings	
	During first 60 hours of sensor wear	After 60 hours of sensor wear	During first 60 hours of sensor wear	After 60 hours of sensor wear
40-80 ^a	62-82%	39%	62-82%	39%
81-120	57-66%	48%	72-84%	66%
a. Agreement to within ± 1.1 mmol/L (20 mg/dL) for glucose readings ≤ 4.4 mmol/L (80 mg/dL)				

Performance of the Guardian RT was evaluated according to the length of time since calibration. This data is not conclusive because of the limited number of data points during the final 3 hours of the 12-hour calibration cycle, i.e., 10. In contrast, 3-hour time bins, earlier in the 12-hour cycle, contained hundreds of data points. This may suggest that calibrations are often required prior to the 12-hour calibration cycle.

Effects of calibration frequency


The average bias when Guardian RT was calibrated ~ 3.5 times a day was -20.5 ± 41 mg/dl (LL: -22.40 mg/dl and UL: -18.63 mg/dl). In those calibrated ~ 5 times a day, the bias was -10.2 mg/dl ± 36 mg/dl (LL: -11.74 mg/dl and UL: -8.66 mg/dl). When comparing Guardian RT units that were calibrated less often to those calibrated more often, the following alarm performance was observed:

- Specificity increased 2-4 percent in the hypoglycemic range and decreased 0-2 percent in the hyperglycemic range
- Sensitivity increased between 5-9 percent across the hyperglycemic range, and decreased 7-16 percent when the alarm was set to 80 mg/dl or below, and decreased 3-7 percent when set between 85 and 100 mg/dl


Stratified error grid analysis also shows better performance in the hypoglycemic range when fewer calibrations are performed, i.e., 62 percent of data points are in Zone A when fewer calibrations were performed, whereas 58 percent were in Zone A when more calibrations were performed.

Time-Elapsed Plots


The plot graph below is a representative example of continuous sensor tracing vs. reference blood glucose reading, where sensor showed excellent performance. The open circles (o) on the graph represent the meter calibration readings. The closed circles (•) represent the reference blood glucose readings, and the solid line (—) represents the sensor glucose value.

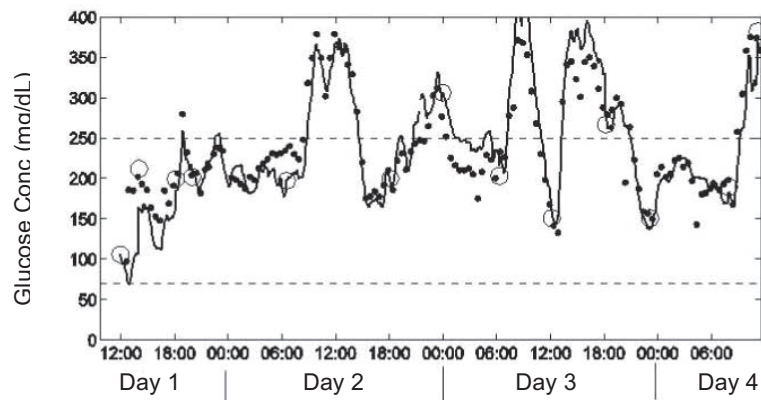

The plot graph below is a representative example of continuous sensor tracing vs. reference blood glucose reading, where sensor showed average (typical) performance. The open circles (o) on the graph represent the meter calibration readings. The closed circles (•) represent the reference blood glucose readings, and the solid line (—) represents the sensor glucose value.

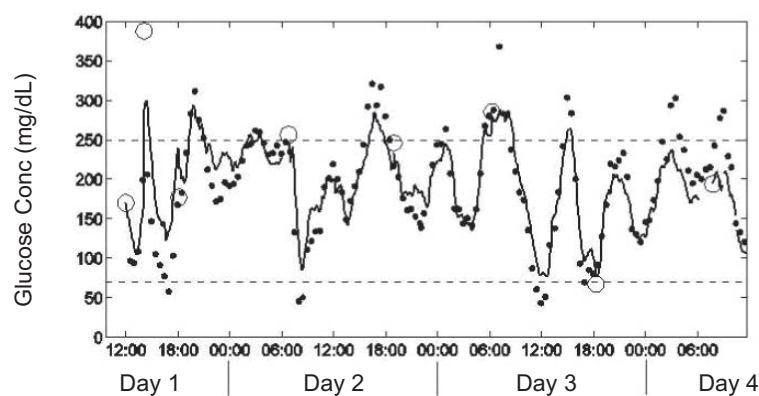

The plot graph below is a representative example of continuous sensor tracing vs. reference blood glucose reading, where sensor showed poor performance. The open circles (o) on the graph represent the meter calibration readings. The closed circles (•) represent the reference blood glucose readings, and the solid line (—) represents the sensor glucose value.

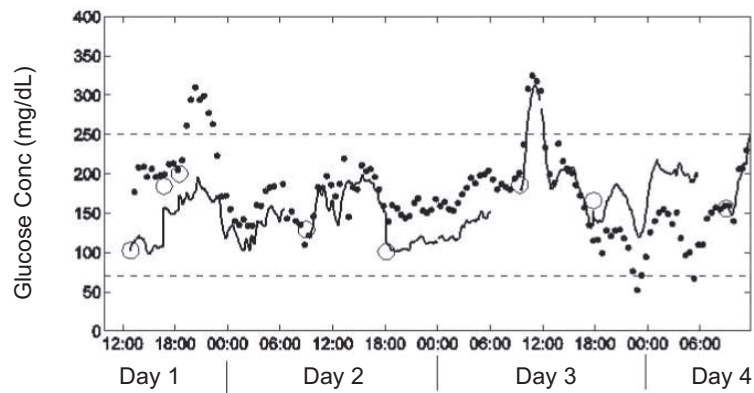

Performance results in children and adolescents

The performance of the Guardian RT in Children and Adolescents was evaluated in a clinical study.¹⁰ Guardian RT results were compared to capillary blood glucose measured by a LifeScan OneTouch® Ultra® Meter. Sixty subjects with Type I diabetes participated in a three-site out-patient study. Subjects ranged in age from 7 to 17 years old. Each subject wore a Guardian RT system, which records a sensor glucose value every five minutes. Subjects were instructed to perform at least seven blood glucose meter measurements each day.

Accuracy of Guardian RT readings

In this study, each blood glucose meter measurement was paired with the corresponding Guardian RT reading (taken every 5 minutes). Pairing was done by selecting the Guardian RT value closest in time to the blood glucose meter result. Agreement was analyzed by comparing paired glucose measurements.

Agreement between the matched pair was estimated by evaluating the difference between the Guardian RT reading and the blood glucose meter measurement. The difference between them was calculated as a percentage of the blood glucose meter measurement (Mean Absolute Percent Difference). The bias was also calculated, and it is defined as the overall difference between the Guardian RT glucose values and the blood glucose meter values. The paired glucose measurements are summarized in the following table

Number of Paired Glucose Measurements	2599
Mean Absolute Percent Difference (± SD)	19.0 ± 19.7%
Bias	-6.0 mg/dL (-0.3 mmol/L)

10Medtronic Diabetes, An Accuracy Evaluation of the Medtronic Diabetes Guardian RT Glucose Monitoring System in Pediatric Subjects with Type 1 Diabetes Mellitus, February, 2006.

The accuracy of the Guardian RT was also evaluated by calculating the percentage of Guardian RT readings within 20% and within 30% of the blood glucose meter reading (or within 20 mg/dL (1.1 mmol/L) in the low glucose range). Results are shown in the following table.

Plasma Glucose Range (mg/dL)	Plasma Glucose Range (mmol/L)	Number of Paired Readings	Percent Within 20%	Percent Within 30%
Overall		2599	68%	82%
40-80*	2.2-4.4	360	51%	----
> 80-120	> 4.4-6.7	482	60%	77%
> 120-240	> 6.7-13.3	1055	74%	90%
> 240	> 13.3	702	75%	89%


*For the Low glucose range, 40-80 mg/dL (2.2-4.4 mmol/L), the value shown is the percent within 20 mg/dL (1.1 mmol/L).

The Clarke Error Grid was used to assess the clinical relevance of the differences between the Guardian RT readings and the comparative blood glucose meter measurements. The Clarke Error Grid divides a correlation plot into 5 zones. See the following table.

Results in zones A and B are considered clinically acceptable, while results in zones C, D, and E are potentially dangerous and, therefore, clinically significant errors. The Clarke Error Grid zones are labeled on the correlation plot.

Zone	Description
A	Clinically accurate, would lead to correct treatment decisions
B	Would lead to benign decisions or no treatment
C	Would lead to over-correction of normal glucose levels
D	Would lead to failure to detect and treat high or low glucose levels
E	Would lead to erroneous treatment decisions

The following figure is a correlation plot of Guardian RT readings versus readings from the reference method, the blood glucose meter. It is overlaid with the Clarke Error Grid. The total number of paired data points is 2599.

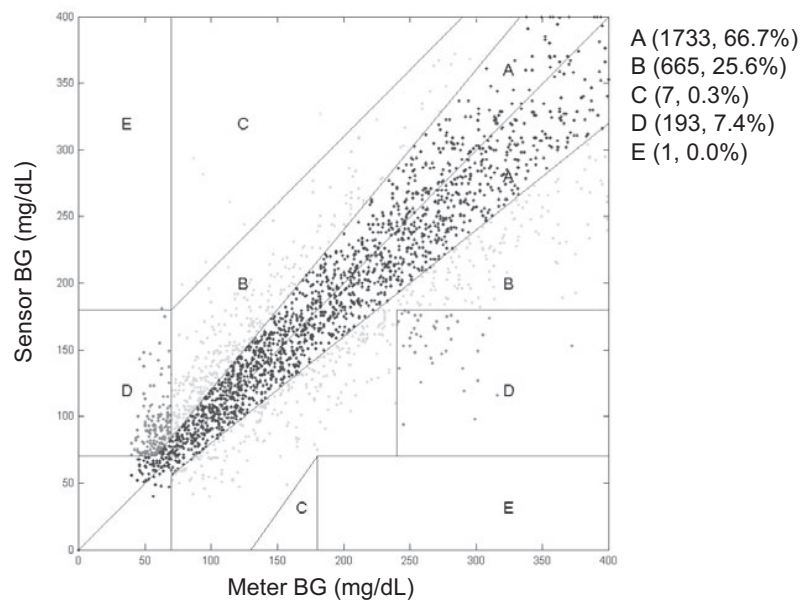

The percent of Guardian RT readings in the above graph are presented in the table below according to the percentage of points falling within each zone (A-E). Results are further broken down (stratified) according to the range of glucose concentrations.

Stratified Clarke Error Grid Analysis

Range of Comparative Glucose Readings (mg/dL)	Total Count	A + B	A	B	C	D	E
40-80	360 (13.9%)	201 (55.8%)	141 (39.2%)	60 (16.7%)	1 (0.3%)	157 (43.6%)	1 (0.3%)
81-120	482 (18.5%)	478 (99.2%)	287 (59.5%)	191 (39.6%)	4 (0.8%)	0 (0%)	0 (0%)
121-240	1055 (40.6%)	1053 (99.8%)	782 (74.1%)	271 (25.7%)	2 (0.2%)	0 (0%)	0 (0%)
240-400	702 (27.0%)	666 (94.9%)	523 (74.5%)	143 (20.4%)	0 (0%)	36 (5.1%)	0 (0%)
Overall	2599 (100.0%)	2398 (92.3%)	1733 (66.7%)	665 (25.6%)	7 (0.3%)	193 (7.4%)	1 (0.0%)

Low and High alerts in children and adolescents

The ability of the Guardian RT to detect high and low glucose levels was measured in the same clinical study.¹¹ Since it is important to set the alert levels in a conservative fashion, the Low Glucose Alert should be set at a value slightly higher than the value of blood glucose you want to detect, and the High Glucose Alert should be set at a value slightly lower than the value of blood glucose you want to detect.

NOTE: Please ask your healthcare professional which low and high alert setting is best for you.

The Low Glucose Alert

The Low Glucose Alert was evaluated for its ability to detect glucose levels at 70 mg/dL (3.9 mmol/L), or below, using the blood glucose meter. As a reference, with the Low Glucose Alert set at 70 mg/dL (3.8 mmol/L), 24% (59/244) of low glucose events were detected by the Guardian RT. Better detection of low blood glucose can be obtained by setting the Low Glucose Alert level higher. For example, setting the Low Glucose Alert at 90 mg/dL (5.0 mmol/L), instead of 70 mg/dL (3.9 mmol/L), increases the ability to detect low blood glucose levels from 24% to 70% (see the following table).

Sometimes the Guardian RT will alert when the blood glucose levels are not low. When the Guardian RT Low Alert was set at 70 mg/dL (3.9 mmol/L) in this study, 48% of the results were considered false alerts (actual blood glucose values are greater than 85 mg/dL (4.7 mmol/L)). This percentage may be exaggerated because blood glucose may be dropping when the Guardian RT alerts.

The table below shows the percent of Low Glucose readings correctly identified by the Guardian RT for specific settings.

Guardian RT Low Alert Setting (mg/dL)	Guardian RT Low Alert Setting (mmol/L)	True Alert Rate ^a	False Alert Rate ^b
70	3.9	24%	48%
80	4.4	52%	46%

a. True Alert Rates are the % of times when the glucose level was at or below the alert setting and the alert sounded.

b. False Alerts Rates are the % of times when the Guardian RT Sensor alarmed but the blood glucose level was greater than the alert setting.

11Medtronic Diabetes, An Accuracy Evaluation of the Medtronic Diabetes Guardian RT Glucose Monitoring System in Pediatric Subjects with Type 1 Diabetes Mellitus, February, 2006.

Guardian RT Low Alert Setting (mg/dL)	Guardian RT Low Alert Setting (mmol/L)	True Alert Rate ^a	False Alert Rate ^b
90	5.0	70%	52%
100	5.6	85%	57%
a. True Alert Rates are the % of times when the glucose level was at or below the alert setting and the alert sounded. b. False Alerts Rates are the % of times when the Guardian RT Sensor alarmed but the blood glucose level was greater than the alert setting.			

Increasing the Low Alert settings will improve the ability to detect low blood glucose events, but it will also increase the frequency of Guardian RT false alerts for blood glucose levels not below the target value. You should consider this trade-off between the improved ability to detect true low blood glucose versus the increased number of false alerts when setting the low alert threshold.

The High Glucose Alert

The High Glucose Alert was evaluated for its ability to detect glucose levels at 250 mg/dL (13.8 mmol/L), or above, using the blood glucose meter. As a reference, with the High Glucose Alert set at 250 mg/dL (13.8 mmol/L), 64% (404/632) of high glucose events were detected by the Guardian RT. Better detection of high blood glucose can be obtained by setting the High Glucose Alert level lower. For example, setting the High Glucose Alert at 190 mg/dL (10.6 mmol/L), instead of 250 mg/dL (13.8 mmol/L), increases the ability to detect high blood glucose levels from 64% to 94% (see the following table).

Sometimes the Guardian RT will alert when the blood glucose levels are not high. When the Guardian RT High Alert was set at 250 mg/dL (13.8 mmol/L) in this study, 13.1% of the results were considered false alerts (actual blood glucose values are less than 225 mg/dL (12.5 mmol/L)). This percentage may be exaggerated because blood glucose may be rising when the Guardian RT alerts. The following table shows the percent of High Glucose readings correctly identified by the Guardian RT for specific settings.

Guardian RT High Alert Setting (mg/dL)	Guardian RT High Alert Setting (mmol/L)	True Alert Rate*	False Alert Rate**
190	10.6	94%	40%
200	11.1	91%	36%

Guardian RT High Alert Setting (mg/dL)	Guardian RT High Alert Setting (mmol/L)	True Alert Rate*	False Alert Rate**
225	12.5	81%	21%
250	13.8	64%	13%

* True Alert Rates are the % of times when the glucose level was at or above the alert setting and the alert sounded.

** False Alerts Rates are the % of times when the Guardian RT Sensor alarmed but the blood glucose level was lower than the alert setting.

Decreasing the High Alert settings will improve the ability to detect high blood glucose events, but it will also increase the frequency of Guardian RT false alerts for blood glucose levels not above the target value. You should consider this trade-off between the improved ability to detect true high blood glucose versus the increased number of false alerts when setting the high alert threshold.

Glossary

A

Adverse reaction - An unexpected, unpleasant or dangerous reaction to a sensor when it is inserted into the body. An adverse reaction may be sudden or may develop over time.

Alarm icon - a solid circle that shows at the top of the screen and the monitor beeps or vibrates periodically until the condition is cleared (see Attention mode).

Alarm Snooze - Once a calibration alert occurs, the monitor will not repeat the alert until after this selected snooze time. This is the setting for the Meter BG Now alert.

Alert icon - an open circle that appears at the top of the screen and the monitor beeps or vibrates periodically to remind you of the condition (see Special mode).

Attention mode - Tells you that an alarm is active or that there is an alarm condition that needs

immediate attention. In the Attention mode, a solid circle shows at the top of the screen and the monitor beeps or vibrates periodically until the condition is cleared.

AUC Limits - The AUC Limits feature allows you to see how long your sensor glucose measurements are outside your target range. Your target range is the area between your Low AUC Limit and your High AUC Limit. When you set the AUC Limits, you select the Low AUC Limit and High AUC Limit you want for your target range.

Autoclave - A chamber for sterilizing with steam under pressure. The original autoclave was essentially a pressure cooker. The steam tightened the lid. The name comes from the Greek *auto*, self + *clavis*, key, or self-locking.

B

BG Units - Blood glucose units used by the monitor (mmol/L or mg/dL).

Blood Glucose - The sugar carried by the blood to provide energy to the human body. The digestive system converts food - protein, fat and carbohydrate - into glucose, which is then absorbed directly into the blood from the small intestine. Also known as dextrose.

C

Calibrate - Check, adjust, or set to a standard (calibrate your monitor).

Cal Reminder - The Cal Reminder reminds the user of an upcoming BG Meter calibration by sending a METER BG BY alert between five minutes and six hours before a new calibration is required (the timing is pre-selected by the user). For example, if the Cal Reminder is set to two hours,

then the METER BG BY alert will occur two hours before the calibration is required.

Circumstance - A condition or fact that has an effect on an action.

Contraindication - A condition that makes a specific treatment or procedure NOT advisable. A contraindication is literally contra- (against) an indication, against something that is indicated as advisable or necessary.

D

Default - A setting or value for your monitor that is assigned automatically by your system. Some default settings cannot be changed; other default settings remain in effect until you change them.

Diabetes mellitus - A group of metabolic diseases featuring high blood glucose (sugar) levels, caused by problems with insulin production. It is generally called Type 1 or Type 2 diabetes and was first identified as a disease about 2,000 years ago. It was called "sweet urine" at that time because high blood glucose levels spill glucose into the urine. Normally, insulin, a hormone produced in the pancreas, controls blood glucose levels, but diabetics either do not produce insulin or their insulin does not work properly. Diabetes cannot

be cured, but it can be controlled through medication, diet and exercise.

Discard - Throw into the trash.

Dressing - A therapeutic or protective material that is put over your sensor.

E

Electromagnetic compatibility - The condition that exists when systems and devices that use electromagnetic energy operate properly without causing or suffering accidental electromagnetic interference to or from other electromagnetic devices.

Electrostatic discharge - The rapid, spontaneous transfer of electrostatic charge induced by a high electrostatic field. The charge usually flows through a spark, known as static discharge, between two objects at different electrostatic levels as they approach each other (for example, when people touch each other after dragging their feet on a carpet).

F

Fingerstick - One method used to test blood glucose by taking a sample of blood for testing from a finger using a lancet or

automatic finger puncture device.

G

Generate - create a result (your monitor will generate a graph).

Glucose - See Blood glucose.

Glucose Limits - These are the high and low Glucose Limits recommended by your healthcare professional. In the Guardian REAL-Time CGM System, a glucose limit is actually a pair consisting of one high limit and one low limit. You can set up the system so that you receive an alert when the sensor glucose measurements reach or exceed your glucose limits.

H

High Glucose alert - This alert occurs when the sensor glucose measurements reach or exceed (go above) your high glucose limit. You have the option of turning this alert on or off.

High Snooze - For the High Glucose, High Predictive, and Rising Rate of change alerts, allows you to set the delay between the first alert and any subsequent alerts. This allows you to avoid getting an alert every five minutes until the condition is corrected.

I

Incineration - The act of burning something completely.

Increment - A small increase or decrease in a measurement.

Indication - A condition that makes a specific medical treatment or procedure advisable or necessary.

Initialize - To set a starting value and/or prepare your monitor, transmitter or sensor for use.

Interruption - When something stops for a while (there was an interruption in communication between your transmitter and your monitor).

Interstitial fluid - The fluid between the cells in the human body. This fluid makes up approximately 16 percent of the weight of the body.

Intravenous - Within or put into a vein; the passing of medicines into a vein through a tube.

Invalidate - An action or a condition that makes a measurement or other information, bad, inaccurate or not usable.

K

Ketone - A chemical created by the human body when there is not enough insulin in the blood.

L

Log - A record of measurements, alarms or other actions by your monitor.

Low Glucose alert - This alert occurs when the sensor glucose measurements reach or exceed (go below) your low glucose limit. You have the option of turning this alert on or off.

Low Snooze - For the Low Glucose, Low Predictive, and Falling Rate of Change alerts, allows you to set the delay between the first alert and any subsequent alerts. This allows you to avoid getting an alert every five minutes until the condition is corrected.

M

Marker - A measurement of a specific event - such as a meal or exercise period, or details about an event - such as the type of insulin you selected the last time you entered information.

Meter - A blood glucose measuring device. If your meter is powered by MWT1 radio frequency (RF) technology, your monitor can be programmed to automatically receive your BG readings from this meter. If your meter is not an RF meter, its data can be entered manually into your monitor.

Missed Data - The monitor will alert if it has not received data from the transmitter for an amount of time that you set.

MWT1 technology - MWT1 is the wireless Radio Frequency (RF) technology that is used to transmit information from the meter to the monitor. You can program your monitor to automatically receive your BG reading from this meter.

O

Occlusive dressing - a bandage that seals a wound from air or bacteria.

P

Predictive alerts - These alerts can be used so that you receive alarms before the sensor glucose measurements reach your glucose limits. This alert predicts when the sensor glucose measurements will reach your glucose limits by using information about your glucose limits and recent sensor glucose measurements.

Preparation - Getting ready to do something.

Prompt - A software message that asks you to do something.

Q

Quaternary ammonium compound - A specific type of antiseptic used to sterilize.

R

Range - The size or area in which something happens.

Rate of Change alerts - These alerts can be used to set up the system so that you receive alerts if the sensor glucose measurements change too rapidly for your needs. You can set the limit for these alerts so that you receive an alert whenever the sensor glucose measurements fall or rise faster than the limit.

Reliable - Always giving the same result or acting the same way.

Reset - To set again, or change the information (reset your glucose values).

Resolve - To find an answer to, or end successfully.

S

Sen-serter - The Sen-serter is an aid for insertion of the Medtronic glucose sensor.

Sensitive - Easily irritated (sensitive skin), or able to measure very slight differences (a sensitive instrument).

Sensor Age - Sensor age is the amount of time, in days and hours, since the sensor was first inserted.

Sensor icon - An icon on your monitor that tells you whether the Sensor feature is On or Off.

Set - To enter or establish a value for your monitor (set your Cal Reminder).

Setting - The position, timing or other way something is set (your High Snooze setting is 15 minutes).

Sharps container - A government-approved container for needles and other "sharp" medical waste containing blood and/or other body fluids. They are usually red in color and made of hard plastic so that needles cannot go through them. They usually have the words "sharps waste" or "biohazard" on them. Check with your healthcare professional about how to properly discard sharps containers.

S/N# - The serial number of the monitor.

Solution - A mixture of two or more substances; often, one is dissolved into another.

Special mode - Tells you that a special feature is active or that a certain alert condition exists. In Special mode, an open circle appears at the top of the screen and the monitor beeps or vibrates periodically to remind

you of the condition. All sensor alerts change from beep/vibration to siren if not cleared within 10 minutes. Special mode does not restrict any monitor functions.

Stabilize - To make steady or keep from changing (stabilize your blood glucose).

Status - The condition of a part of your system (battery status, alarm/alert status).

Sterile - Free from substances or germs that could cause infection.

Sterile gauze - A thin medical bandage with loose weave that is free from substances or germs that could cause infection.

Submerge - To place under water or cover completely with water.

T

Timeout - If you do not enter information or do something in a specific length of time, your monitor automatically stops what it is doing and goes to the HOME screen.

Transmitter - The device that sends a signal to your monitor with information from your sensor.

Transmtr Batt - The status of the transmitter battery. Possible values are "Good," "Low," and "Bad."

Transmtr ID - The identification number of the transmitter that you are using.

Transmtr Ver - The software version of the transmitter that you are using.

Index

Numerics

3-hour graph 61
6-hour graph 62
12-hour graph 63
24-hour graph 64

A

A (alarm) condition 94
accessories
 CareLink USB 4
 holster 4
 leather case 4
 meter 4
 monitor clip 4
ADD ALARM screen 79
Alarm Clock 78
 default settings 112
Alarm Clock reminders
 deleting 79
 reviewing 79
alarm conditions 19
Alarm conditions
 A (alarm) 94
 Batt out limit 94
 Button error 94
 Check settings 94
 E (Error) 94
 Failed batt test 95

Low battery 92
Off no power 95
Reset 95
Weak Battery 95
alarm history
 specifications 110
ALARM HISTORY screen 95, 101
alarm messages and codes
 specifications 110
ALARM OPTION screen 79
Alarm Silence
 default settings 112
 options 70
ALARM SILENCE screen 70
Alarm Snooze
 default settings 112
 setting 35
Alarm/alert icons 18
Alarm/alert types 19
alarms
 basic types 89, 90, 91, 92
 Cal error 90
 Change transmtr 90
 Falling Rate of Change 92
 High glucose 91
 High Predictive 91
 Lost sensor 89
 Low glucose 91
 Low Predictive 91
 Low transmtr 89
 Meter BG by 91
 Meter BG now 91

Rising Rate of Change 92
Sensor end 90
Sensor error 92
Weak signal 89
Change sensor 90
conditions 89, 94
E- (error) 94
 how to respond 88
 reviewing recent alarms 95
Alert Type
 default settings 112
alert types
 battery power usage for Beep
 versus Vibrate 26
 Beep 26
 specifications 110
 Vibrate 26
alerts
 Cal Error 90
 for alert conditions 92
 low battery alert 92
Attention mode 20
AUC limits 35
AUC Limits
 default settings 113
 entering 35

B

backlight 15, 18, 92, 102
 specifications 110

- Batt Out Limit alarm 12
- Batt out limit alarm condition 94
- battery 18, 19, 27, 101, 102, 103, 105
 - AAA Energizer alkaline 12
 - cap 13
 - compartment 103, 106
 - icon 17, 18
 - icons 17
 - installing 12
 - monitor 22, 67, 95, 101, 105
 - rechargeable or carbon zinc
 - battery warning 12
 - recommendations 12
 - transmitter 69
- battery:
 - monitor 102
- Battery icon 17, 18
- BG units 72
- BG Units
 - default settings 113
- bleeding 45
- Block
 - default settings 113
- Block feature 80
- BLOCK OPTION screen 80
- blood glucose
 - entering incorrect number 90
 - rapid change 90
- blood glucose
 - measurements 35, 54, 55
 - non-calibration meter,
 - examples of 55
- blood glucose units
 - setting 28
- blood glucose units (BG Units)
 - default settings 113
- Button error alarm condition 94
- buttons 11, 17, 21, 102
 - ESC button 22
 - using sharp objects with 11

C

- Cal error alarm condition 90
- Cal Error alert 90
- Cal Reminder
 - default settings 113
- CAL REMINDER screen 28, 36
- Calibrating your system 51
- calibrating your system
 - Capture Event feature 54
 - entering exercise
 - information 57
 - information types that must
 - be manually logged into the
 - system 54
- CAPTURE EVENT 20
- Capturing Events 54
- CARB UNITS screen 56
- CareLink Personal 85
- CareLink USB 4
- Change sensor alarm condition 90
- Change transmtr alarm
 - condition 90
- charger 97
 - connecting with transmitter 41
- charging
 - transmitter 41
- Check settings alarm condition 94
- Checking battery status 97
- Checking the sensor status 68
- Checking the status of the
 - system 67
- Checking the transmitter status 69
- choosing a sensor insertion site 43
- choosing an insertion site
 - ones to avoid 43
- Cleaning the monitor 106
- Cleaning the Sen-serter 107
- Cleaning the system 106
- Cleaning the transmitter 107, 108
- clear settings 83

- clear the alarm 88
- clearing the alarm 20
- Connecting the transmitter to
 - the sensor 46
- connecting transmitter to
 - charger 41
- continuous graph display 37
- contraindications 7
- Conventions used in this guide 5
- current settings
 - reviewing 11, 72

D

- Daily totals 75
- daily totals
 - specifications 112
- Damage due to radiation 8
- default screen
 - specifications 112
- default settings
 - Alarm Clock 112
 - Alarm Silence 112
 - Alarm Snooze 112
 - Alert Type 112
 - AUC Limits 113
 - BG Units 113
 - Block 113
 - blood glucose units (BG
 - Units) 113
 - Cal Reminder 113
 - Glucose Alarms 114
 - Glucose Limits 114
 - Glucose Limits 1 114
 - Glucose Limits 2 through 8 115
 - High Glucose Limit 115
 - High Snooze 115
 - Language 116
 - Lock Keypad 116
 - Low Glucose Limit 115
 - Low Snooze 116

- Meter Option 116
- Missed Data 117
- Predictive Alarms 117
- sensor 117
- sensor glucose graph
 - timeout 118
- specifications 112
- Time/Date setup 118
- DELETE ALARM screen 80
- deleting Alarm Clock reminders 79
- disconnecting
 - transmitter from sensor 72
- disposing
 - introducer needle 45

E

- E (Error) alarm condition 94
- EDIT SETTINGS
 - screen 27, 32, 35, 37
- electromagnetic emissions 121
- Electromagnetic environment 126
- electromagnetic immunity 122
- emergency kit 2
- emissions 121
- ENTER BG screen 48, 53, 55
- ENTER FOOD screen 57
- ENTER INSULIN screen 56
- ENTER MEAL screen 56
- entering exercise information 57
- entering first meter blood glucose
 - procedure 47
- Entering information about
 - carbohydrates 55
- Entering insulin injection
 - information 55
- Entering non-calibration meter
 - BG measurements 55
- Examples of real-time sensor
 - glucose graphs 61
- EXERCISE DURATION screen 57

EXERCISE INTENSITY screen 57

F

- Failed Batt Test alarm 12
- Failed batt test alarm condition 95
- Falling Rate of Change alarm
 - condition 92
- Find lost sensor 101

G

- Glucose Alarms
 - default settings 114
 - turning on 28
- GLUCOSE ALARMS screen 28
- Glucose Limits
 - default settings 114
- Glucose Limits 1
 - default settings 114
- Glucose Limits 2 through 8
 - default settings 115
- glucose monitoring settings
 - entering 27, 28, 30, 32, 34, 35, 36, 37, 38
 - reviewing your settings 38
 - setting AUC limits 35
 - setting Glucose Limits 28
 - setting Glucose Limits,
 - guidelines 30
 - setting Glucose Limits,
 - setting start time 30
 - setting Glucose Limits,
 - switching limits off 28, 30
 - setting Predictive glucose
 - alarm 32
 - setting Rate of Change
 - alarms 34

- setting sensor glucose
 - graph timeout 37
- setting the Cal Reminder 36
- turning on Glucose
 - Alarms 28

- graphs 60
 - differentiating between
 - current and earlier events
 - in 59
 - opening 59, 60
 - real-time glucose 58
 - real-time sensor glucose
 - levels 60
 - selecting 15
- Guardian REAL-Time CGM System
 - effect of airport security
 - systems on 8
 - indications 1
 - introduction 1
 - items included with 3
 - monitor 3
 - transmitter 3

H

- High glucose alarm condition 91
- High Glucose Limit
 - default settings 115
 - setting 29
 - specifications 111
- High Predictive alarm 33
- High Predictive alarm condition 91
- High Snooze
 - default settings 115
 - setting 31
- High Snooze alarm
 - setting 32
- history
 - sensor updates 71
- How to make glucose alarms
 - silent 70

I

- Icons on the monitor
 - Alert 18
 - Battery 17, 18
- injection insulin types 55
- inserting sensor
 - Sen-serter 44
 - stopping bleeding 45
- inserting the sensor
 - sites to avoid 43
- insertion sites
 - adverse reactions at 42
 - choosing 43
 - preparing 43
- insulin 55
- INSULIN TYPE screen 55
- Intermediate Acting insulin 56
- introducer needle 44
 - disposing 45
 - removing 45

L

- Language
 - default settings 116
- Lock keypad 80
- Lock Keypad
 - default settings 116
- LOCK KEYPAD screen 80
- locking the keypad 80
- Long Acting insulin 56
- Lost sensor alarm condition 89
- low battery 18
 - affects your alert type 27
- Low battery alarm condition 92
- Low glucose alarm condition 91
- Low Glucose Limit
 - default settings 115
 - setting 29

- Low Glucose
 - Limit:specifications 111
- Low Predictive alarm condition 91
- Low Snooze
 - default settings 116
- Low Snooze alarm
 - setting 32
- Low Snooze feature 32
- Low transmtr alarm condition 89

M

- MAIN MENU 15, 20
- measurements
 - to calibrate system 51, 71
 - system calibration, guidelines 51
- menus 20
- Meter BG by alarm condition 91
- Meter BG Now alarm 35
- Meter BG now alarm condition 91
- meter IDs
 - add, delete, review 50
- meter option 48
 - rules 49
- Meter Option
 - default settings 116
- METER OPTION screen 50
- METER STATS screen 77, 78
- MiniLink charger 40
- MiniLink transmitter 97
- MiniLink transmitter
 - specifications 127
- Missed Data
 - default settings 117
- Missed Data feature 89
 - description 37
- Missed Data option 37
- Mixed 50/50 insulin 56
- Mixed 70/30 insulin 56

- modes
 - Attention 94
 - Normal 19
 - Special 19
- monitor 11
 - current settings 11
 - icons 17, 18
 - Alarm 18
 - Battery 17, 18
 - Sensor 17
 - Time of day 18
 - menus 20
 - CAPTURE EVENT 20
 - MAIN MENU 20
 - SENSOR MENU 20
 - operating modes 19, 20
 - Attention 20
 - Normal 19
 - Special 19
 - removing and storing 105
 - screens 17, 22
 - STATUS 22
 - uses for 17
 - scroll bar 18
 - storing 105
- monitor,
 - using in water 103
- monitor battery status
 - Low 67
 - Normal 67
 - Off 67
- Monitor data management
 - feature 75

N

- needle guard 44
- NEW SENSOR screen 47
- Next Cal 68
- Normal mode 19

O

Off no power alarm condition 95
on an aircraft
 entering meter blood glucose
 measurements 9
OTHER screen 58

P

Precautions 7
 extreme temperatures 7
 sensor 8
Predictive Alarms 117
PREDICTIVE ALARMS LOW/HIGH
 screen 33
preparing sensor insertion sites 44

R

Rate of Change alarms
 sensitivity and alarm
 frequency 34
reading the real-time glucose
 graphs 58
real-time glucose graphs
 reading 58
real-time sensor glucose levels 60
rechargeable or carbon zinc
 battery warning 12
removing
 introducer needle 45
 needle guard 44
 sensor 73
removing the transmitter and
 sensor 73
Reset alarm condition 95
Restore settings 83
Restore Settings feature 82

REVIEW ALARMS screen 79
REVIEW SETTINGS screen 38, 72
reviewing Alarm Clock
 reminders 79
reviewing current settings 72
reviewing recent alarms 95
RF communications 126
RF features
 meter 48
Rising Rate of Change alarm
 condition 92

S

Safety Warnings 7
Save Settings feature 82
screens
 ADD ALARM 79
 ALARM HISTORY 95
 ALARM OPTION 79
 BLOCK OPTION 80
 DELETE ALARM 80
 LOCK KEYPAD 80
 METER STATS 77, 78
 REVIEW ALARMS 79
 SENSOR STATS 77, 78
 SETTINGS HISTORY 84
Selftest 81, 102, 103
Sen-serter 42
 inserting sensor 44
 removing from sensor 45
Sen-serter, storing the 105
sensor
 alarms 90
 Change sensor 90
 Change transmr 90
 best sites to use for
 insertion 43
 checking status 68
 choosing an insertion site 43
 default settings 117

 description 42
 disconnecting from
 transmitter 72
 inserting 45
 bleeding 45
 insertion sites to avoid 43
 removing 73
 stabilizing 90
 Transmitter ID 36, 37
 viewing alarm history 93
Sensor
 alarms 88, 89, 90, 91, 92
 Cal error 90
 conditions 89
 Falling Rate of Change 92
 High glucose 91
 High Predictive 91
 how to respond 88
 Lost sensor 89
 Low glucose 91
 Low Predictive 91
 Low transmr 89
 Meter BG by 91
 Meter BG now 91
 Rising Rate of Change 92
 Sensor end 90
 Sensor error 92
 Weak signal 89
 alarms/alerts 87
 understanding 87
 attaching the transmitter 40
 procedure 40
 connecting to the
 transmitter 47
 starting 47
 update history, viewing
 calibration values 71
Sensor Age 68
Sensor Alarm History 93
Sensor alerts 89
sensor alerts
 Cal Error 90

- Sensor Demo 38
- Sensor end alarm condition 90
- Sensor error alarm condition 92
- Sensor function 20
- sensor glucose graph timeout
 - default settings 118
- sensor glucose graphs
 - continuous display 37
- sensor glucose levels
 - graphs 60
- sensor glucose measurements 87
- Sensor icon 17
- Sensor Isig 68
- SENSOR MENU 20
- SENSOR ON/OFF screen 27, 38
- SENSOR READY 2 HRS screen 47
- SENSOR STATS screen 77, 78
- sensor status 68
 - Next Cal 68
 - Sensor Age 68
 - Sensor Isig 68
 - Transmtr Batt 69
 - Transmtr ID 69
 - Transmtr Ver 69
- SENSOR STATUS screen 68, 69
- SENSOR UPDATE HISTORY screen 71
- sensor updates
 - viewing history 71
- SET ALARM SNOOZE screen 35
- SET AUC LIMITS screen 35
- SET CAL REMINDER screen 36
- SET FALL RATE LIMIT screen 34
- SET GLUCOSE LIMITS 1 screen 30
- SET GLUCOSE LIMITS 3 screen 31
- SET GRAPH TIMEOUT screen 37
- SET HIGH SNOOZE screen 32
- SET HOUR screen 24
- SET LOW SNOOZE screen 32
- SET MISSED DATA screen 37
- SET RISE RATE LIMIT screen 34
- SET START TIME 2 screen 31
- SET TIME SENSITIVITY LOW/HIGH screen 33
- SET TRANSMITTER ID screen 36
- setting
 - High Snooze 31
 - Low Snooze alarm 32
 - time and date 24
- setting the Alarm Snooze 35
- Setting the sensor glucose graph timeout 37
- setting up your sensor and transmitter
 - basic steps 39, 40, 47
 - attaching the transmitter to your body 40
 - connecting the sensor to the transmitter 47
 - entering first meter blood glucose 47
 - starting the sensor 47
- settings
 - entering 24, 26, 27, 28, 30, 32, 34, 35, 36, 37, 38
 - basic steps 24
 - glucose monitoring settings 27
 - glucose monitoring settings, reviewing your settings 38
 - glucose monitoring settings, setting AUC limits 35
 - glucose monitoring settings, setting Glucose Limits, guidelines 30
 - glucose monitoring settings, setting Glucose Limits, setting start time 30
 - glucose monitoring settings, setting
- Predictive glucose alarm 32
- glucose monitoring
 - settings, setting Rate of Change alarms 34
- glucose monitoring
 - settings, setting sensor glucose graph timeout 37
- glucose monitoring
 - settings, setting the Cal Reminder 36
- glucose monitoring
 - settings, turning on Glucose Alarms 28
- selecting alert type 26
- setting time and date 24
- entering, basic steps 23
- glucose monitoring
 - 27, 28, 30, 32, 34, 35, 36, 37, 38
 - reviewing your settings 38
 - setting AUC limits 35
 - setting Glucose Limits 28
 - setting Glucose Limits, guidelines 30
 - setting Glucose Limits, setting start time 30
 - setting Predictive glucose alarm 32
 - setting Rate of Change alarms 34
 - setting sensor glucose graph timeout 37
 - setting the Cal Reminder 36
 - turning on Glucose Alarms 28
 - reviewing current 72
- SETTINGS HISTORY screen 84
- Short Acting insulin 55
- Special mode 19
- Starting the sensor 47

- starting the sensor
 - procedure 47
- status
 - sensor 68
- STATUS screen 101
- stopping bleeding 45
- storage temperature, monitor 128
- Storing the MiniLink
 - transmitter 105
- storing the monitor 105
- Storing the Sen-serter 105
- Storing the system 105
- supplies, ordering 4
- System alarm conditions 94
- system alarm/alert history
 - viewing 95
- system settings
 - clearing 82, 83
 - restoring 83
 - saving 82

T

- tester 99
- Testing the transmitter 99
- the HOME screen 17
- time and date
 - setting 24
- Time of day icon 18
- TIME/DATE SET screen 24
- Time/Date setup
 - default settings 118
- TIME/DATE SETUP screen 24
- Transmitter
 - attaching to your body 40
 - connecting the sensor to 47
- transmitter
 - before attaching 45
 - charging 41
 - charging complete 42
 - checking status 69

- connecting with charger 41
- description 42
- disconnecting from sensor 72
- lights 41
- Transmitter ID
 - setting 36, 37
- transmitter status 69
- transmitter 69
- transmitter, storing the
 - MiniLink 105
- transmitter, testing the 99
- Transmtr Batt 69
- Transmtr ID 69
- Transmtr Ver 69

U

- unlocking the keypad 81
- User safety
 - contraindications 7
 - Damage due to radiation 8
 - indications 6
 - Precautions 7
 - warnings 7
- User settings 82
- User Settings screen 103
- User Settings, viewing your 84
- using the Capture Event
 - feature 54
- Using your system in water 73
- UTILITIES MENU 21

V

- Viewing sensor alert history 93
- viewing sensor update history 71
- Viewing your current settings 72
- viewing your system alarm/alert
 - history 95

- viewing your User Settings
 - history 84

W

- Weak Battery alarm 12
- Weak Battery alarm condition 95
- Weak signal alarm condition 89
- What to do when you get an
 - alarm/alert 88
- Where to insert the sensor 43



6025254-012C 012907