



ARTISAN[®]
TECHNOLOGY GROUP

Your **definitive** source
for quality pre-owned
equipment.

Artisan Technology Group

(217) 352-9330 | sales@artisanTG.com | artisanTG.com

Full-service, independent repair center

with experienced engineers and technicians on staff.

We buy your excess, underutilized, and idle equipment

along with credit for buybacks and trade-ins.

Custom engineering

so your equipment works exactly as you specify.

- Critical and expedited services
- In stock / Ready-to-ship
- Leasing / Rentals / Demos
- ITAR-certified secure asset solutions

Expert team | Trust guarantee | 100% satisfaction

All trademarks, brand names, and brands appearing herein are the property of their respective owners.

1

Find the **GE Healthcare / Amersham Biosciences MegaBACE 4500** at our website: **Click *HERE***

1029

Alamar Exhibit 1029
Alamar v. Olink

MegaBACE 4000

Instrument Administrator's Guide
Version 3.2



MegaBACE, Molecular Dynamics, and SNUPe are trademarks of Amersham Biosciences Limited.

Amersham and Amersham Biosciences are trademarks of Amersham plc.

Cimarron is a trademark of Cimarron Software, Inc.

Microsoft and Windows are trademarks of Microsoft Corporation.

The MegaBACE DNA Analysis System is for research purposes only. It is not intended or approved for diagnosis of disease in humans or animals.

All goods and services are sold subject to the terms and conditions of sale of the company within the Amersham Biosciences group that supplies them. A copy of these terms and conditions is available on request.

© Amersham Biosciences Corp 2002—All rights reserved.



December 2002

Notice to purchaser: limited license

The MegaBACE instrument is a confocal scanning system licensed under US Patent Numbers 5,091,652 and 5,274,240, and corresponding foreign patents and patent applications, including any continuations, continuations-in-part, and subdivisions and the like.

The instrument is also an Authorized DNA Sequencer. It is authorized under one or more US Patent Numbers 4,849,513; 5,171,534; 5,015,733; 5,118,800; 5,161,507; 5,118,802; 4,855,225; and 5,366,860, and corresponding foreign patents and patent applications. The purchase of this instrument includes limited, non-exclusive rights under the subject patents to use this instrument for sequencing and fragment length analysis when used with Authorized Reagents. The use of this instrument with Authorized Reagents provides a limited license to perform DNA sequencing and fragment length analysis in accordance with the label rights accompanying such reagents. Purchase of this instrument does not itself convey to the purchaser a complete license to perform DNA sequencing and fragment length analysis under the subject patents. Authorized reagents may be obtained from licensed vendors, or reagents may be authorized under separate license arrangements from Applied Biosystems. The above patent rights are granted solely for research and other uses that are not unlawful. No other licenses are granted expressly, impliedly, or by estoppel.

Further information on purchasing licenses to perform DNA sequencing and fragment length analysis may be obtained by contacting the Director of Licensing at Applied Biosystems, 850 Lincoln Center Drive, Foster City, California 94404.

Applied Biosystems does not guarantee the performance of this instrument.

Amersham Biosciences is a licensed vendor for authorized reagents.

Amersham Biosciences UK Limited Amersham Place Little Chalfont
Buckinghamshire England HP7 9NA

Amersham Biosciences AB SE-751 84 Uppsala Sweden

Amersham Biosciences Corp 800 Centennial Avenue PO Box 1327 Piscataway
NJ 08855 USA

Amersham Biosciences Europe GmbH Munzinger Strasse 9 D-79111 Freiburg
Germany

Amersham Biosciences (SV) Corp 928 East Arques Avenue Sunnyvale
CA 94085-4520 USA

Table of contents

Preface

About this guide xiii

Related publications. xiii

Assumptions xiii

Assistance xiv

Part one General system configuration

Chapter 1 Introduction to the MegaBACE 4000 system software

1.1 The MegaBACE 4000 system software components 1-1

1.2 Overview of the Instrument Control Manager features 1-2

1.3 Overview of the sequencing setup 1-2

1.3.1 Administrator tasks for sequencing 1-3

1.3.2 The system initialization files for sequencing 1-3

1.3.3 About data collection with automatic base calling 1-4

1.3.4 The parameter template files for sequencing. 1-4

1.3.5 The plate setup data files (.psd) for sequencing 1-5

1.3.6 Workflow alternatives for the sequencing plate setup. 1-5

1.3.7 How the software uses the sequencing parameters 1-7

1.4 Overview of the STR genotyping setup 1-8

1.4.1 Administrator tasks for STR genotyping. 1-8

1.4.2 The system initialization files for STR genotyping. 1-9

1.4.3 The parameter template files for STR genotyping. 1-9

1.4.4 The plate setup data files (.psd) for STR genotyping 1-10

1.4.5 Workflow alternatives for the STR genotyping plate setup. . . . 1-10

1.4.6 How the software uses the STR genotyping parameters. 1-12

1.5 Overview of the SNP genotyping setup 1-13

1.5.1 Administrator tasks for SNP genotyping. 1-13

1.5.2 The system initialization files for SNP genotyping. 1-14

1.5.3 The parameter templates for SNP genotyping 1-14

Table of contents

1.5.4	The plate setup data files (.psd) for SNP genotyping.	1-15
1.5.5	Workflow alternatives for the SNP genotyping plate setup . . .	1-15
1.5.6	SNP genotyping plate sets and plate attributes	1-17
1.5.7	How the software uses the SNP genotyping parameters	1-18
1.6	File format overview.	1-20
1.7	File storage overview	1-20

Chapter 2 Specifying basic system preferences

2.1	Using the MegaBACE.ini file	2-1
2.1.1	Modifying the MegaBACE.ini file.	2-3
2.1.2	General field names.	2-4
2.2	Specifying a default application	2-5
2.3	Enabling or disabling the edit mode	2-6
2.4	Specifying how plate setup data (.psd) files are used	2-7
2.4.1	Specifying a new .psd file folder location	2-8
2.4.2	Specifying that .psd files are required.	2-8
2.5	Using automated bar coding	2-9
2.5.1	Specifying the plate ID length.	2-9
2.5.2	Specifying the plate ID terminator.	2-10
2.6	Suppressing file creation for empty or bad capillaries.	2-10

Chapter 3 Specifying the chemistry parameters

3.1	How the software uses the Chemistry.ini file.	3-1
3.2	About specifying the chemistry	3-1
3.3	The fields in the Chemistry.ini file.	3-2
3.4	Adding a new set of chemistry parameters.	3-3

Chapter 4 Creating templates for the instrument control parameters

4.1	About instrument control templates	4-1
4.2	Using the Instrument Control window to create a template	4-4
4.3	Using a text file to create an instrument control template (.icp)	4-5

Part two Setting up the system for sequencing

Chapter 5 Specifying preferences for sequencing runs

5.1	The sequencing fields in the MegaBACE.ini file	5-1
-----	--	-----

5.2	Setting default parameter templates	5-6
5.2.1	Setting the default plate setup template.	5-6
5.2.2	Setting the default instrument control template	5-6
5.3	Configuring the automatic base calling feature	5-7
5.4	Turning on or off the automatic base calling feature	5-7
5.5	Setting up the automatic file export feature.	5-9
5.5.1	Specifying the file formats for export	5-9
5.5.2	The FASTA file format	5-10
5.5.3	Changing the SCF version.	5-10
5.5.4	Using subfolders for the exported files.	5-10
5.5.5	Using the file export options	5-11
5.6	Settings in the Basecall.ini file	5-11
5.7	Specifying which base caller to use	5-12
5.7.1	Specifying a default base caller.	5-12
5.7.2	About specifying a base caller(s) in the plate definition	5-13

Chapter 6 Creating templates for sequencing plate setup parameters

6.1	About creating the plate definitions.	6-1
6.1.1	Methods for creating the plate definitions	6-1
6.1.2	The methods for creating plate definitions in the Plate Setup window	6-2
6.2	About plate setup templates for sequencing	6-3
6.3	Using the Plate Setup window to create a plate setup template	6-5
6.4	Creating a plate setup template in a text file	6-6
6.5	Using a template to suppress .rsd file creation for empty capillaries	6-7
6.6	Setting and clearing the default plate setup template	6-8
6.6.1	Using the Templates menu to set a default template	6-9
6.6.2	Clearing the default plate setup template.	6-9

Chapter 7 Creating plate setup data files (.psd) for sequencing

7.1	Using plate setup data files (.psd) for sequencing.	7-1
7.1.1	About the plate ID.	7-1
7.1.2	Using a plate ID that is different from the scanned ID	7-3

Table of contents

7.2	Creating a plate setup data file (.psd)	7-4
7.2.1	The .psd field names for sequencing	7-5
7.2.2	The .psd file value constraints	7-6
7.3	Specifying a base caller in a .psd file	7-7
7.3.1	Designating a base caller for the plate	7-7
7.3.2	Specifying a base caller for a selected well on a plate	7-8
7.4	Annotating well-specific information in a .psd file	7-8
7.5	Designating the sample names and file names in the .psd file	7-9
7.6	Using a .psd file to suppress .rsd file creation for empty capillaries . .	7-10

Part three Setting up the system for STR genotyping

Chapter 8 Specifying preferences for STR genotyping

8.1	The STR genotyping fields in the MegaBACE.ini file	8-1
8.2	Setting STR genotyping as the default application	8-3
8.3	Setting default parameter templates for STR genotyping	8-3
8.3.1	Setting the default plate setup template for STR genotyping . .	8-4
8.3.2	Setting the default instrument control template for STR genotyping	8-4

Chapter 9 Creating templates for STR genotyping plate setup
parameters

9.1	About creating the plate definitions	9-1
9.1.1	Methods for creating the plate definitions	9-1
9.1.2	The methods for creating plate definitions in the Plate Setup window	9-2
9.2	About plate setup templates for STR genotyping	9-3
9.3	Using the Plate Setup window to create a template	9-5
9.4	Creating a plate setup template in a text file	9-6
9.5	Using a template to suppress .rsd file creation for empty capillaries . .	9-6
9.6	Setting and clearing the default plate setup template	9-8
9.6.1	Using the Templates menu to set a default template	9-8
9.6.2	Clearing the default plate setup template	9-8

Chapter 10 Creating plate setup data files (.psd) for STR genotyping

10.1	Using plate setup data files (.psd) for STR genotyping	10-1
10.1.1	About the plate ID	10-1
10.1.2	Using a plate ID that is different from the scanned ID	10-2
10.2	Creating a plate setup data file (.psd)	10-4
10.2.1	The .psd field names.	10-5
10.2.2	The .psd file value constraints	10-6
10.3	Annotating well-specific information in a .psd file	10-6
10.4	Designating the sample names and file names in the .psd file.	10-7
10.5	Using a .psd file to suppress .rsd file creation for empty capillaries . .	10-8

Part four Setting up the system for SNP genotyping

Chapter 11 Specifying preferences for SNP genotyping

11.1	The SNP genotyping fields in the MegaBACE.ini file	11-1
11.2	Setting SNP genotyping as the default application.	11-5
11.3	Setting default parameter templates for SNP genotyping	11-6
11.3.1	Setting the default plate set setup template for SNP genotyping	11-6
11.3.2	Setting the default instrument control template for SNP genotyping	11-7
11.4	Specifying limits for the injections.	11-7
11.4.1	Specifying a maximum number of injections	11-7
11.4.2	Specifying the injection timeout period	11-8
11.5	Enforcing the order in which an operator injects the plates	11-9
11.6	Prompting the operator for plate-specific information	11-9
11.6.1	Providing missing SNP marker names	11-10
11.6.2	Providing missing plate IDs	11-12
11.6.3	Copying sample names from a previous injection.	11-12

Chapter 12 Creating templates for SNP genotyping plate set setup parameters

12.1	About creating the plate set definitions	12-1
12.1.1	Methods for creating the plate set definitions	12-2
12.1.2	The methods for creating plate set definitions in the Plate Set Setup window	12-2

Table of contents

12.2	About the plate set setup parameters for SNP genotyping	12-4
12.3	Using the Plate Set Setup window to create a template	12-5
12.4	Creating a plate set setup template in a text file (.tpl)	12-6
12.5	Using a template to suppress .rsd file creation for empty capillaries .	12-7
12.6	Setting and clearing the default plate set setup template	12-9
12.6.1	Using the Templates menu to set a default template	12-9
12.6.2	Clearing the default plate set setup template.	12-9

Chapter 13 Creating plate setup data files (.psd) for SNP genotyping

13.1	Creating plate setup data files (.psd) for SNP genotyping	13-1
13.1.1	About the types of .psd files	13-1
13.1.2	Creating a plate setup data file (.psd)	13-4
13.2	Specifying the plate set definition in a master .psd file	13-7
13.2.1	Specifying the plate set ID	13-7
13.2.2	Defining the list of plates in the plate set.	13-8
13.2.3	Defining the number of injections	13-9
13.3	Using a .psd file to suppress .rsd file creation for empty capillaries. .	13-9
13.4	Specifying the attributes for a plate in a .psd file	13-10
13.4.1	Specifying the plate ID	13-10
13.4.2	Specifying the SNP marker name	13-10
13.4.3	Specifying sample names for a given plate	13-11

Part five Protocol modification and file storage

Chapter 14 Bypassing or overriding a protocol

14.1	Bypassing the selected protocol.	14-2
14.2	Overriding the contents of a protocol	14-3

Chapter 15 File storage

15.1	Raw sample data storage.	15-1
15.1.1	Raw sample data file content	15-1
15.1.2	Naming and storage conventions (sequencing and STR genotyping only)	15-1
15.1.3	Naming and storage conventions (SNP genotyping only). . .	15-2

15.2	Base-called sample data storage (sequencing only)	15-3
15.2.1	Naming and storage conventions for the base-called and exported files	15-3
15.2.2	Naming and storage conventions for the base-called run folders	15-5
15.3	System initialization file storage	15-6
15.4	Parameter configuration file storage	15-7

Part six **Appendixes**

Appendix A Header formats

A.1	The standard header format	A-1
A.2	Viewing and editing the data in the extended header	A-3
A.3	The basic extended header format	A-4
A.4	The extended header labels in sequencing files	A-4
A.4.1	The extended header in the sequencing .rsd files	A-4
A.4.2	The bar code	A-6
A.4.3	The extended header in base-called sample data files	A-6
A.5	The extended header labels in STR genotyping files	A-7
A.6	The extended header labels in SNP genotyping files	A-7
A.6.1	Additional extended header labels for SNP genotyping	A-7

Appendix B Protocol for self-describing data

B.1	Introduction	B-1
B.2	The 1-byte header	B-1
B.3	The length value	B-2
B.4	The data buffer	B-2
B.5	Notes about the protocol	B-3

Appendix C Launching a cleanup script after a run

C.1	Using the CLEANUP COMMAND	C-1
C.2	Batch file examples	C-3
C.2.1	Displaying the full set of arguments after the run	C-3
C.2.2	Coping the raw sample data files to another folder.	C-5

Glossary

Preface

About this guide

The *MegaBACE 4000 Instrument Administrator's Guide* describes how to configure the MegaBACE™ system software for your laboratory. This guide describes how to use the system initialization files, how to create plate setup and instrument control templates, and how to create the plate setup data files (.psd) to automate the workflow for the instrument operator.

If you have instrument operator responsibilities for the MegaBACE 4000 system in addition to being an administrator, see the *MegaBACE 4000 Instrument User's Guide* for additional guidelines.

Related publications

In addition to the *MegaBACE 4000 Instrument Administrator's Guide*, the following publications are available for the MegaBACE 4000 system:

- *MegaBACE 4000 Instrument User's Guide* describes how to use the MegaBACE system to perform runs and how to use the instrument protocols. This guide also provides instructions on maintaining the instrument and guidelines on troubleshooting problems.
- *MegaBACE 4000 Site Preparation Guide* provides instructions for setting up the installation site for the MegaBACE 4000 instrument. Reading the site preparation guide is a prerequisite for the installation of the MegaBACE 4000 system.
- The MegaBACE analysis software user's guides and the Help available within the software provide instructions on how to use the software to analyze the data collected from the MegaBACE 4000 instrument.

Assumptions

The software-related instructions in this user's guide assume you have basic computer skills. You should be familiar with the Microsoft™ Windows™ graphical user interface. If you do not have these skills, refer to the documentation or the Help for Windows.

Assistance

When calling for assistance, be prepared to supply the serial number of your instrument. The serial number is located on the lower right side of the MegaBACE instrument. For contact by phone or fax, please use one of the numbers below.

Asia Pacific

Tel: +852 2811 8693
Fax: +852 2811 5251

Australasia

Tel: +61 2 9899 0999
Fax: +61 2 9899 7511

Austria

Tel: 01 576 0616 10
Fax: 01 576 0616 27

Belgium

Tel: 0800 73 888
Fax: 03 272 1637

Canada

Tel: +1 800 463 5800
Fax: +1 800 567 1008

Central, East, and Southeast Europe

Tel: +43 1 982 3826
Fax: +43 1 985 8327

Denmark

Tel: 45 16 2400
Fax: 45 16 2424

Finland & Baltics

Tel: +358 (0)9 512 39 40
Fax: +358 (0)9 512 17 10

France

Tel: 01 69 35 67 00
Fax: 01 69 41 96 77

Germany

Tel: 0761 4903 291
Fax: 0761 4903 405

Italy

Tel: 02 27322 1
Fax: 02 27302 212

Japan

Tel: +81 3 5331 9336
Fax: +81 3 5331 9370

Latin America

Tel: +55 11 3933 7300
Fax: +55 11 3933 7315

Middle East and Africa

Tel: +30 (1) 96 00 687
Fax: +30 (1) 96 00 693

Netherlands

Tel: 0165 580 410
Fax: 0165 580 401

Norway

Tel: 2318 5800
Fax: 2318 6800

Portugal

Tel: 21 417 70 35
Fax: 21 417 31 84

Russia & other C.I.S. & N.I.S.

Tel: +7 (095) 232 0250, 956 1137
Fax: +7 (095) 230 6377

Southeast Asia

Tel: +60 3 8024 2080
Fax: +60 3 8024 2090

Spain

Tel: 93 594 49 50
Fax: 93 594 49 55

Sweden

Tel: 018 612 1900
Fax: 018 612 1910

Switzerland

Tel: 01 802 81 50
Fax: 01 802 81 51

UK

Tel: 0800 616928
Fax: 0800 616927

USA

Tel: +1 800 526 3593
Fax: +1 877 295 8102

Web site

<http://www.amershambiosciences.com>

Part one

General system configuration

Chapter 1 Introduction to the MegaBACE 4000 system software

The MegaBACE 4000 DNA Analysis System is a high-throughput automated gene analysis system with 384-capillary capacity. The system consists of reagents for sample preparation, the MegaBACE instrument, and the software.

This chapter provides an overview of the system software. The topics are—

- The MegaBACE 4000 system software components (section 1.1)
- Overview of the Instrument Control Manager features (section 1.2)
- Overview of the sequencing setup (section 1.3)
- Overview of the STR genotyping setup (section 1.4)
- Overview of the SNP genotyping setup (section 1.5)
- File format overview (section 1.6)
- File storage overview (section 1.7)

1.1 The MegaBACE 4000 system software components

The MegaBACE 4000 instrument control software consists of the—

- **Instrument Control Manager**—Contains the graphical interface that allows you to control the instrument.
- **Host Scan Controller**—Provides the communication between the MegaBACE 4000 instrument and the Instrument Control Manager.

Analysis software is available for use with the MegaBACE 4000 instrument to support various types of applications, including:

- Sequencing analysis
- Microsatellite genotyping, also known as short-tandem-repeat (STR) genotyping, and fragment analysis
- Single-nucleotide-polymorphism (SNP) genotyping

The MegaBACE Header Editor allows you to view and edit the data in the extended header of the raw sample data files (.rsd) (section A.2).

1.2 Overview of the Instrument Control Manager features

The Instrument Control Manager allows you to—

- Specify the plate definitions in plate setup data files (.psd) that the operator can import at the start of the Inject Samples and Run protocol. Alternatively, you or the operator can use the Plate Setup window to create plate definitions. The operator can then select these precreated plate definitions at the start of the Inject Samples and Run protocol.
- Specify the default plate setup parameters and instrument parameters using templates. The MegaBACE 4000 software includes templates (default) for the recommended parameters.
- Enable or disable the edit mode. You can allow the instrument operator to edit or prevent the ability to edit the plate setup and instrument parameters for running a plate.
- Change applications. The software supports sequencing, STR genotyping (or fragment analysis), and SNP genotyping.
- (SNP genotyping only) Serially inject multiple sample plates into the same set of capillaries for a single run.
- Collect data.
- (Sequencing only) Automatically call bases on the collected data and export the data to other file formats after base calling.
- Change the data storage location, and store the raw run folders in application-specific locations (sequencing, STR genotyping, and SNP genotyping).

1.3 Overview of the sequencing setup

This section provides an overview of the administrator's tasks and the system setup for the sequencing application.

1.3.1 Administrator tasks for sequencing

Table 1-1 provides an administrator’s task overview for the sequencing application on the MegaBACE 4000 system. The administrator tasks vary depending on whether your laboratory uses plate setup data files (.psd). The .psd files automate the operator’s workflow for providing the plate definition for a run.

Table 1-1. Administrator’s task overview for sequencing

Task	Frequency	Reference
To tailor the Instrument Control Manager to the laboratory’s workflow, you can modify the— • MegaBACE.ini file • Basecall.ini file • Chemistry.ini file	Once, when the laboratory first receives the system, or whenever the system needs to be reconfigured	Chapters 2, 3, and 5
The software includes default template files for the— • Plate setup parameters (.tpl files) • Instrument control parameters (.icp files) You can create additional templates based on the needs of your laboratory, and you can specify defaults.	Once, when the laboratory first receives the system, or whenever you change the run conditions for multiple plates	Chapters 4 and 6
(Optional) You can create plate setup data files (.psd) to automate how an operator provides the plate definition for a run. For each .psd file you create, you provide the operator with a .psd file name or bar code to import the file.	For each plate you run	Chapter 7

1.3.2 The system initialization files for sequencing

To meet the needs of your laboratory, you might need to modify the system initialization files. The system software includes three initialization files that you can modify—

- **MegaBACE.ini file**—Determines the settings the Instrument Control Manager software uses to control the operations of the MegaBACE 4000 instrument (chapters 2 and 5). For example, the MegaBACE.ini file can specify—
 - The default application (sequencing, STR genotyping, or SNP genotyping) and the default plate setup and instrument control templates for each application
 - Whether an operator can manually edit the plate setup or instrument control parameters

- Whether the software performs base calling automatically after a run is completed
- **Chemistry.ini file**—Determines the chemistry parameter sets available in the Instrument Control Manager software. See chapter 3 for details on how to use the Chemistry.ini file.
- **Basecall.ini file**—Determines the base callers available in the Instrument Control Manager software, and specifies a default base caller. See chapter 5 for details on how to use the Basecall.ini file.

1.3.3 About data collection with automatic base calling

The Instrument Control Manager allows you to perform continuous base calling and file export after data collection without operator intervention. To do this, you specify automatic base calling, select a default base caller, and specify the file formats for the base-called data (SCF, FASTA, ASCII text, or ABD). For a description of how to set up the base calling feature, see chapter 5.

1.3.4 The parameter template files for sequencing

The parameter templates are optional files that can be used to define the plate setup parameters and instrument control parameters that an operator can import at run time. The software includes templates that are ready for you to use or modify. You can also create new templates. The Instrument Control Manager uses two kinds of parameter templates—

- **Plate setup template file (.tpl)**—You can use a .tpl file for multiple plates to specify only the plate setup parameters that appear in the Plate Setup window of the Instrument Control Manager. The software includes a selection of templates for sequencing. See chapter 6 for a description of how to create a plate setup template.
- **Instrument control template file (.icp)**—You can use an .icp file to specify the instrument control parameters for multiple runs. The instrument control parameters define the instrument run conditions and the matrix fill and flush cycles. The software includes a default template for sequencing (Normal.icp). See chapter 4 for a description of how to create an instrument control template.

Important**1.3.5 The plate setup data files (.psd) for sequencing**

For the software to import the .psd file, the text or bar code that an operator enters must match the file name (case insensitive) of the .psd file.

An operator can automatically import the plate definition from a .psd file. A .psd file can specify the plate ID, plate setup parameters, and other plate- and well-level attributes, such as sample names and file names.

The .psd files can be used, either alone or in combination with a .tpl file, to automate the operator's task of defining a plate setup. If you use a .psd file to specify all the plate setup parameters, you do not need to use a .tpl file. Similarly, if you use a .tpl file to specify all the plate setup parameters for multiple plates, you do not need to use a .psd file unless you want to specify additional parameters for a particular plate.

The parameters entered last, whether entered manually in the Plate Setup window or automatically using a .psd or .tpl file, overwrite any previously entered parameters. The parameters that you specify in the .psd file become part of the extended header of each raw sample data (.rsd) file created during the run. Section 7.1 provides examples of .psd files for sequencing.

1.3.6 Workflow alternatives for the sequencing plate setup

During a run on the MegaBACE 4000 instrument, the Instrument Control Manager prompts the operator for the appropriate steps. The software includes protocols for rinsing the capillaries, filling the capillaries with sieving matrix, performing a prerun, and performing the sample injection for the electrophoresis run.

For each run, the operator must provide a plate definition. A plate definition consists of the—

- **Plate ID**—The name of the plate. The software uses the plate ID to name the raw run folder that stores the data from the run.
- **Plate setup parameters**—The electrophoresis parameters, the chemistry parameters, the sample names, plate comments, and other parameters.

The workflow that the operator follows depends on how you configure the Instrument Control Manager and whether your laboratory uses plate setup data files (.psd). Figure 1-1 shows the basic workflows, which are:

- **Automated workflow**—For the most automated workflow, the operator can import the plate definition from a .psd file at the start of the Inject Samples and Run protocol. As the administrator, you create the .psd file (chapter 7) and provide the operator with the file name or the bar code to enter for the run.
- **Alternative workflow**—
 - You or the operator can use the Plate Setup window to create the plate definition before the run. The values for the plate setup parameters can be provided from a template, imported from a .psd file, or manually entered in the Plate Setup window. See chapter 6 of this guide for a description of how to create a template. See the *MegaBACE 4000 Instrument User's Guide* for a description of how to use the Plate Setup window to create plate definitions.
 - At the start of the Inject Samples and Run protocol, the operator selects the plate ID for the precreated plate definition.

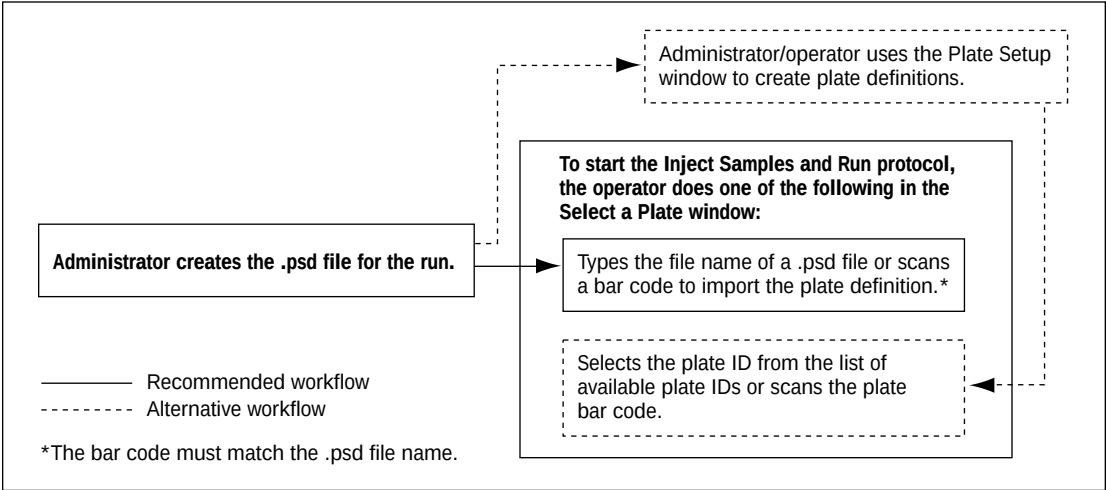


Figure 1-1. Workflow alternatives for providing the sequencing plate definition for a run.

1.3.7 How the software uses the sequencing parameters

See figure 1-2 for an overview of how the Instrument Control Manager uses the parameters specified in the .psd files, in the parameter templates (.tpl and .icp), and by keyboard entry.

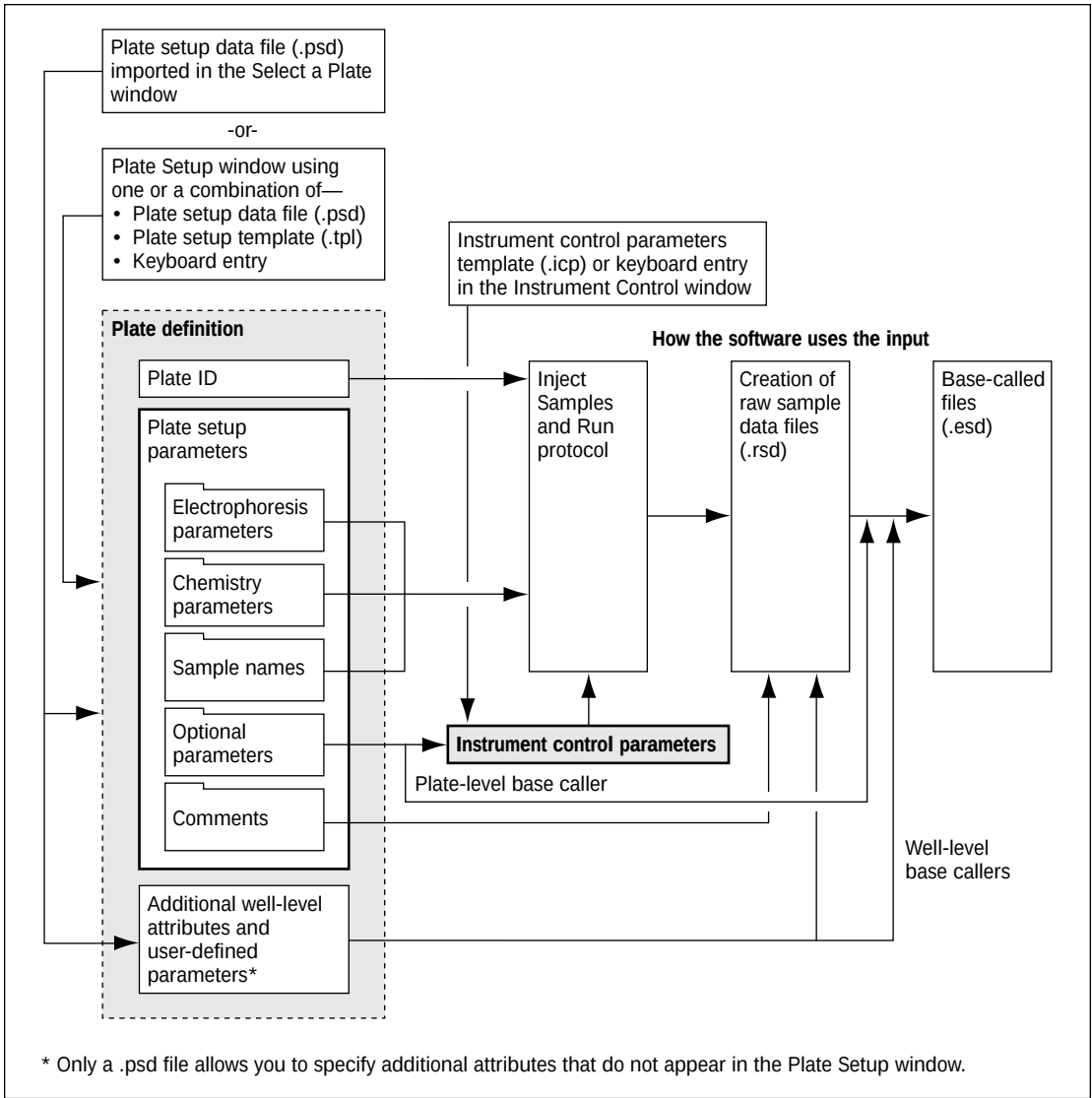


Figure 1-2. Overview of how the software uses the parameters specified for a sequencing run.

1.4 Overview of the STR genotyping setup

This section provides an overview of the workflow and the system setup for the STR genotyping application.

1.4.1 Administrator tasks for STR genotyping

Table 1-2 provides an administrator's task overview for the STR genotyping application on the MegaBACE 4000 system. The administrator tasks vary depending on whether your laboratory uses plate setup data files (.psd). The .psd files automate the operator's workflow for providing the plate definition for a run.

Table 1-2. Administrator's task overview for STR genotyping

Task	Frequency	Reference
To tailor the Instrument Control Manager to the laboratory's workflow, you can modify the— • MegaBACE.ini file • Chemistry.ini file	Once, when the laboratory first receives the system, or whenever the system needs to be reconfigured	Chapters 2, 3, and 8
The software includes default template files for the— • Plate setup parameters (.tpl files) • Instrument control parameters (.icp files) You can create additional templates based on the needs of your laboratory, and you can specify defaults.	Once, when the laboratory first receives the system, or whenever you change the run conditions for a set of plates	Chapters 4 and 9
(Optional) You can create plate setup data files (.psd) to automate how an operator provides the plate definition. For each .psd file you create, you provide the operator with a .psd file name or bar code to import the file.	For each plate you run	Chapter 10

1.4.2 The system initialization files for STR genotyping

To meet the needs of your laboratory, you might need to modify the system initialization files. The system software includes two initialization files that you can modify—

- **MegaBACE.ini file**—Determines the settings the Instrument Control Manager uses to control the operations of the MegaBACE 4000 instrument (chapters 2 and 8). For example, the MegaBACE.ini file can specify—
 - The default application (sequencing, STR genotyping, or SNP genotyping) and the default plate setup and instrument control templates for each application
 - Whether an operator can manually edit the plate setup or instrument control parameters
- **Chemistry.ini file**—Determines the chemistry parameter sets available in the Instrument Control Manager software. See chapter 3 for details on how to use the Chemistry.ini file.

1.4.3 The parameter template files for STR genotyping

The parameter templates are optional files that can be used to define the plate setup parameters and instrument control parameters that an operator can import at run time. The software includes templates that are ready for you to use or modify, or you can create new ones. The Instrument Control Manager uses two kinds of parameter templates—

- **Plate setup template file (.tpl)**—You can use a .tpl file for multiple plates to specify only the plate setup parameters that appear in the Plate Setup window of the Instrument Control Manager. The software includes a default template for genotyping (StdGenotyping.tpl). See chapter 9 for a description of how to create a plate setup template for STR genotyping.
- **Instrument control parameter template file (.icp)**—You can use an .icp file to specify the instrument control parameters for multiple runs. The instrument control parameters define the instrument run conditions and the matrix fill and flush cycles. The software includes a default template for genotyping (Genotyping.icp). See chapter 4 for a description of how to create an instrument control template.

Important

1.4.4 The plate setup data files (.psd) for STR genotyping

For the software to import the .psd file, the text or bar code that an operator enters must match the file name (case insensitive) of the .psd file.

An operator can automatically import the plate definition from a .psd file. A .psd file can specify the plate ID, plate setup parameters, and other plate- and well-level attributes, such as sample names and file names.

The .psd files can be used, either alone or in combination with a .tpl file, to automate the operator's task of defining a plate. If you use a .psd file to specify all the plate setup parameters, you do not need to use a .tpl file. Similarly, if you use a .tpl file to specify all the plate setup parameters for multiple plates, you do not need to use a .psd file unless you want to specify additional parameters for a particular plate.

The parameters entered last, whether entered manually in the Plate Setup window or automatically using a .psd or .tpl file, overwrite any previously entered parameters. The parameters that you specify in the .psd file become part of the extended header of each raw sample data (.rsd) file created during the run. Section 10.1 provides examples of .psd files for STR genotyping.

1.4.5 Workflow alternatives for the STR genotyping plate setup

During a run on the MegaBACE 4000 instrument, the Instrument Control Manager prompts the operator for the appropriate steps. The software includes protocols for rinsing the capillaries, filling the capillaries with sieving matrix, performing a prerun, and performing the sample injection for the electrophoresis run.

For each run, the operator must provide a plate definition. A plate definition consists of the—

- **Plate ID**—The name of the plate. The software uses the plate ID to name the raw run folder that stores the data from the run.
- **Plate setup parameters**—The electrophoresis parameters, the chemistry parameters, the sample names, plate comments, and other parameters.

The workflow that the operator follows depends on how you configure the Instrument Control Manager and whether your laboratory uses plate setup data files (.psd). Figure 1-3 shows the workflows, which are—

- **Automated workflow**—For the most automated workflow, the operator can import the plate definition from a .psd file at the start of the Inject Samples and Run protocol. As the administrator, you create the .psd file (chapter 10) and provide the operator with the file name or the bar code to enter for the run.
- **Alternative workflow**—
 - You or the operator can use the Plate Setup window to create the plate definition before the run. The values for the plate setup parameters can be provided from a template, imported from a .psd file, or manually entered in the Plate Setup window. See chapter 9 of this guide for a description of how to create a template. See the *MegaBACE 4000 Instrument User's Guide* for a description of how to use the Plate Setup window to create plate definitions.
 - At the start of the Inject Samples and Run protocol, the operator selects the plate ID of the precreated plate definition.

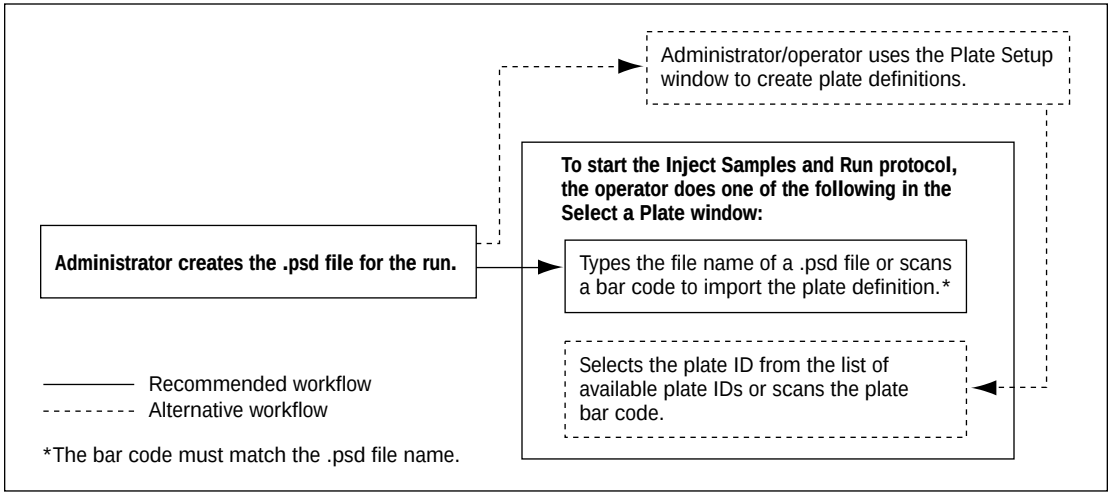


Figure 1-3. Typical STR genotyping workflows for providing the plate definition for a run.

1.4.6 How the software uses the STR genotyping parameters

See figure 1-4 for an overview of how the Instrument Control Manager uses the parameters specified in the .psd files, in the parameter templates (.tpl and .icp), and by keyboard entry.

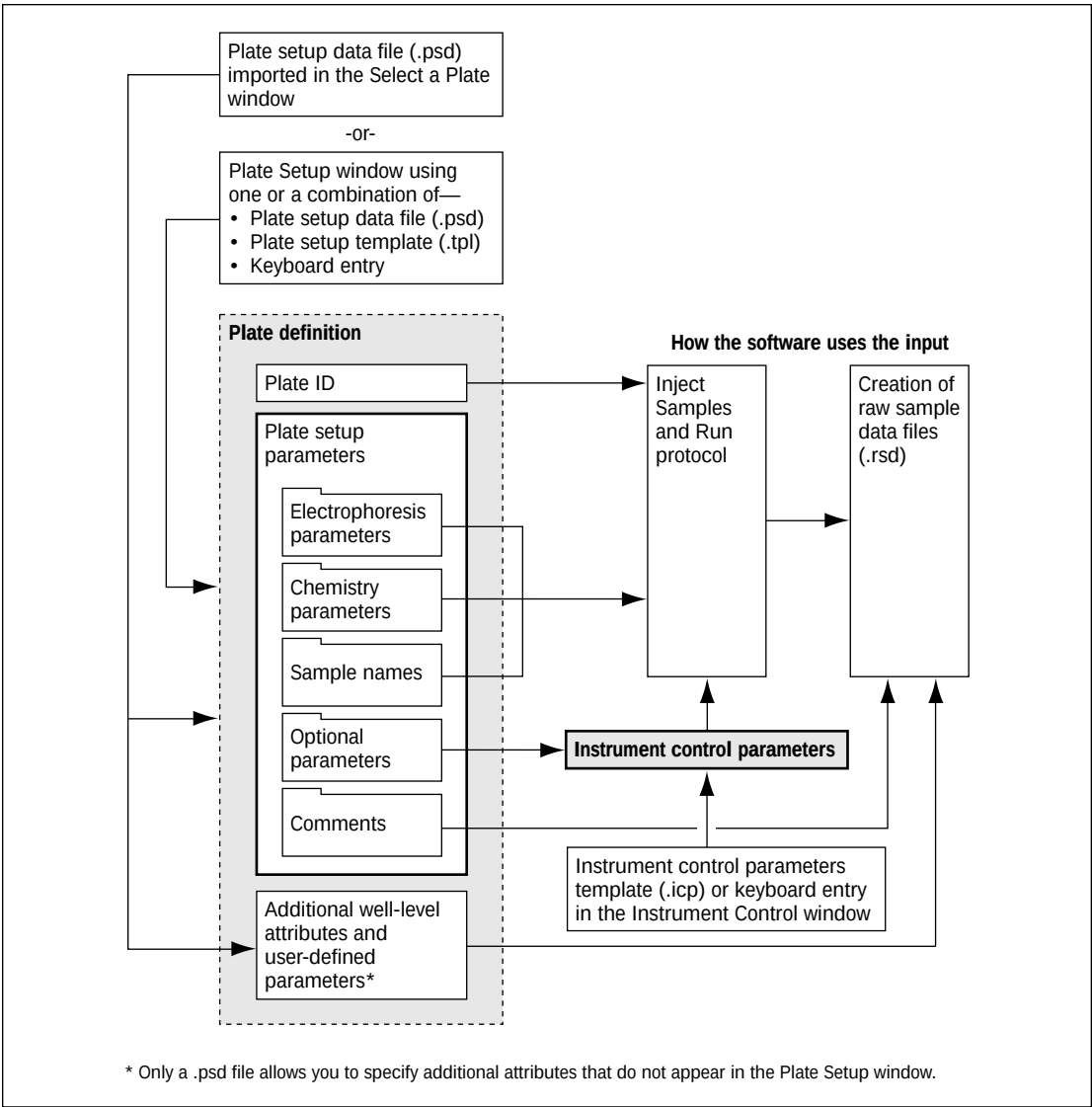


Figure 1-4. Overview of how the software uses the parameters specified for an STR genotyping run.

1.5 Overview of the SNP genotyping setup

This section provides an overview of the administrator’s tasks and the system setup for the SNP genotyping application.

1.5.1 Administrator tasks for SNP genotyping

Table 1-3 provides an administrator’s task overview for the SNP genotyping application on the MegaBACE 4000 system. The administrator tasks vary depending on whether your laboratory uses plate setup data files (.psd). To automate the workflow during a multi-injection run, .psd files are recommended.

Table 1-3. Administrator’s task overview for SNP genotyping

Task	Frequency	Reference
To tailor the Instrument Control Manager to the laboratory’s workflow, you can modify the— - MegaBACE.ini file - Chemistry.ini file	Once, when the laboratory first receives the system, or whenever the system needs to be reconfigured	Chapters 2, 3, and 11
The software includes default template files for the— - Plate set setup parameters (.tpl files) - Instrument control parameters (.icp files) You can create additional templates based on the needs of your laboratory, and you can set defaults.	Once, when the laboratory first receives the system, or whenever you change the run conditions for multiple plates	Chapters 4 and 12
(Optional) You can create plate setup data files (.psd) to automate how an operator provides the plate set* definition and the attributes for each plate injected during a run. For each .psd file you create, you provide the operator with a .psd file name or bar code to import the file.	For each plate set and the associated plates in the set	Chapter 13

* A plate set includes all the plates injected during a multi-injection run for the SNP genotyping application. See section 1.5.6 for an overview of the plate set definition and the plate set setup parameters.

1.5.2 The system initialization files for SNP genotyping

To tailor the SNP genotyping application for your laboratory, you can modify the system initialization files. The system software includes two initialization files that you can modify—

- **MegaBACE.ini file**—Determines the settings the Instrument Control Manager uses to control the operation of the MegaBACE instrument.

For example, the MegaBACE.ini file can specify—

- The default application (sequencing, genotyping, or SNP genotyping) and the default plate setup and instrument control templates for each application
- Whether an operator can manually edit the plate set setup or instrument control parameters
- Whether the software enforces the injection order of the plates in a plate set
- Whether the software prompts the operator for a plate ID and SNP marker name, or whether the software provides the option of copying the sample names from the previous injection

- **Chemistry.ini file**—Determines the chemistry parameter sets available in the Instrument Control Manager software. The Chemistry.ini file includes a chemistry parameter set for the SNuPe™ chemistry.

1.5.3 The parameter templates for SNP genotyping

The Instrument Control Manager allows you to use templates to specify two different kinds of parameters—

- **Plate set setup template (.tpl)**—You can use a .tpl file for multiple plate sets to specify only the parameters that appear in the Plate Set Setup window of the Instrument Control Manager. The plate set setup parameters are required for the plate set definition (section 1.5.6). The software includes a default plate set setup template for SNP genotyping (StdSNP.tpl). See chapter 12 for a description of how to create a plate set setup template for SNP genotyping.
- **Instrument control template (.icp)**—You can use an .icp file to specify the instrument control parameters for multiple runs. The instrument control parameters define the instrument run conditions and the matrix fill and flush cycles. The software includes a default instrument control template for SNP genotyping (SNP_Typing.icp). See chapter 4 for a description of how to create a default instrument control template.

The parameter templates can automate the operator tasks for the SNP genotyping application.

Important

1.5.4 The plate setup data files (.psd) for SNP genotyping

For the software to import the .psd file, the text or bar code that an operator enters must match the file name (case insensitive) of the .psd file.

For each multi-injection run, you can provide two types of .psd files (chapter 13)—

- **Master .psd file**—Allows you to specify the plate set definition, including a list of the plates to be included in the plate set. You also use the master .psd file to specify the plate-specific attributes for the first plate to be injected.
- **Plate-specific .psd files**—Allow you to specify the plate-specific attributes for each injection.

For example, for a multi-injection run consisting of nine plates (injections), you would prepare—

- **One master .psd file**—To define the plate set and provide the attributes for the first plate to be injected.
- **Eight additional plate .psd files**—To provide the attributes for the second injection and all the subsequent injections in the run.

1.5.5 Workflow alternatives for the SNP genotyping plate setup

The SNP samples are prepared using the SNuPe genotyping kit and the SNuPe multi-injection marker kit. During a single run, the operator can serially inject multiple SNuPe sample plates into the same set of capillaries (multi-injection run). The Instrument Control Manager prompts the operator for the appropriate steps, including rinsing the capillaries, filling the capillaries with sieving matrix, performing a prerun, and performing the multiple sample injections for the electrophoresis run.

For each multi-injection run, the operator must provide—

- A plate set definition, which includes the plate set ID and the setup parameters for the plate set
- The plate-specific attributes for each injection

The workflow that the operator follows depends on how you configure the Instrument Control Manager and whether your laboratory uses plate setup data files (.psd). Figure 1-5 shows the workflow alternatives.

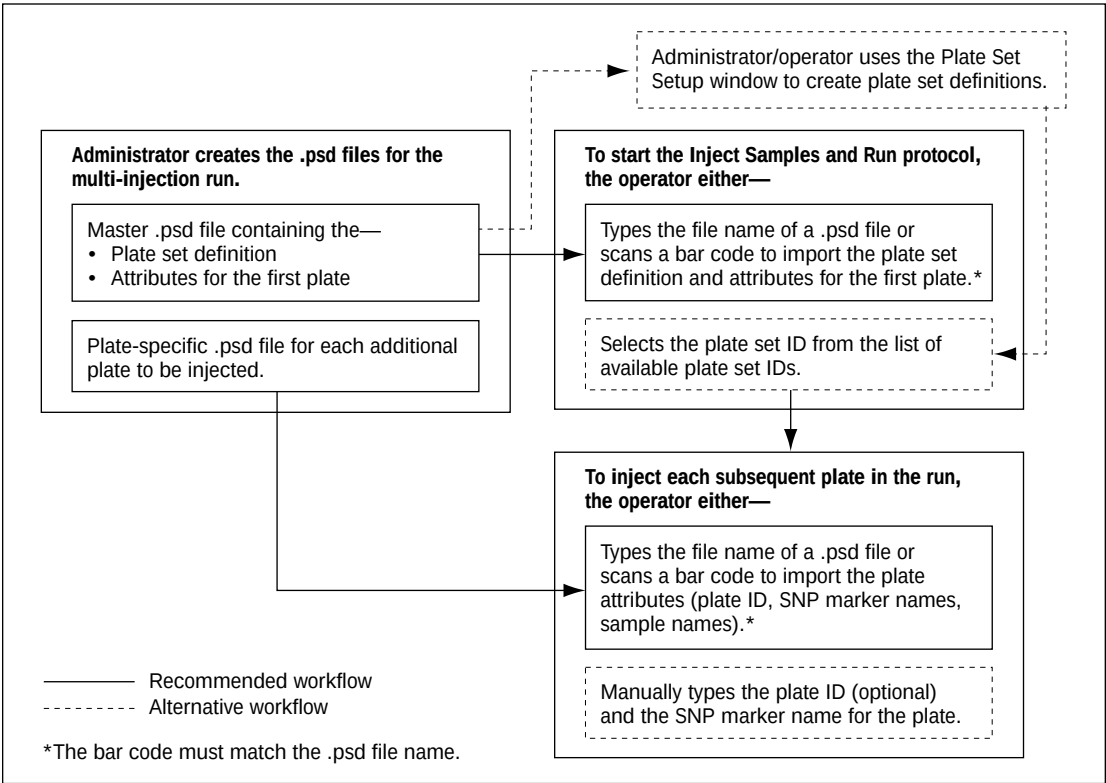


Figure 1-5. Workflow alternatives for providing the plate set definition and the plate-specific attributes for SNP genotyping runs.

As figure 1-5 shows, the workflow alternatives include—

- **Automated multi-injection workflow**—For the most automated workflow, the operator can import the required information from .psd files during the Inject Samples and Run protocol. As the administrator, you create the .psd files and provide the operator with the file names or bar codes to enter during the run:
 - A master .psd file containing the plate set definition and attributes for the first plate
 - Additional .psd files to import the plate-specific attributes for each subsequent sample injection
- **Alternative multi-injection workflow**—
 - You or the operator can use the Plate Set Setup window to create the plate set definition before the run. The values for the plate set setup parameters can be provided from a template, imported from a .psd file, or manually entered in the window. See the *MegaBACE 4000 Instrument User's Guide* for a description of how to use the Plate Set Setup window.
 - At the start of the Inject Samples and Run protocol, the operator selects the plate set ID for the precreated plate set definition.
 - For each injection, the operator can manually enter a few plate-specific attributes when the plate is injected.

1.5.6 SNP genotyping plate sets and plate attributes

The plate set definition for SNP genotyping

A plate set definition consists of the—

- **Plate set ID**—The name for the set of plates injected in a multi-injection run. The software uses the plate set ID to name the raw run folder that stores the data from the run. **Note:** If the plate set ID is missing, the software uses the ID of the plate that is injected first as the ID for the entire plate set. If the plate ID is also missing, the software uses the master .psd file name.
- **Plate set setup parameters**—The electrophoresis parameters, the chemistry parameters, the file names, plate set comments, and other parameters shared by all the plates injected in a single run. The values for these parameters can be imported from a master .psd file, selected from a plate set setup template, or manually entered in the Plate Set Setup window.

- **(Optional) List of plates**—The list of bar codes representing the plates contained in the plate set. To specify the list of plates and the injection order for a plate set, you must use a master .psd file.

The plate-specific attributes for SNP genotyping

The plate-specific attributes for the multi-injection run can include—

- **(Optional) Plate ID**—You can provide a plate ID in the .psd file that the operator imports for the injection or you can allow the operator to manually type the plate ID for the injection.
- **SNP marker names**—The Instrument Control Manager requires at least one SNP marker name for each injection. To specify different SNP marker names for different wells on a plate, the software requires a .psd file for the injection. Alternatively, you can allow the operator to manually type one SNP marker name for each injected plate.
- **(Optional) Sample names**—To specify the sample names for a plate, the Instrument Control Manager requires a .psd file for the sample plate injection.
- **(Optional) Plate comments**—To specify plate comments and other user-defined plate attributes, the software requires you to use a .psd file for the injection. **Note:** Comments about the plate set are included as part of the plate set setup parameters.

If a .psd file is not used for a given injection, the operator can specify only the plate ID and one SNP marker name for the injection.

1.5.7 How the software uses the SNP genotyping parameters

See figure 1-6 for an overview of how the Instrument Control Manager uses the parameters specified in the .psd files, in the parameter templates (.tpl and .icp), and by keyboard entry.

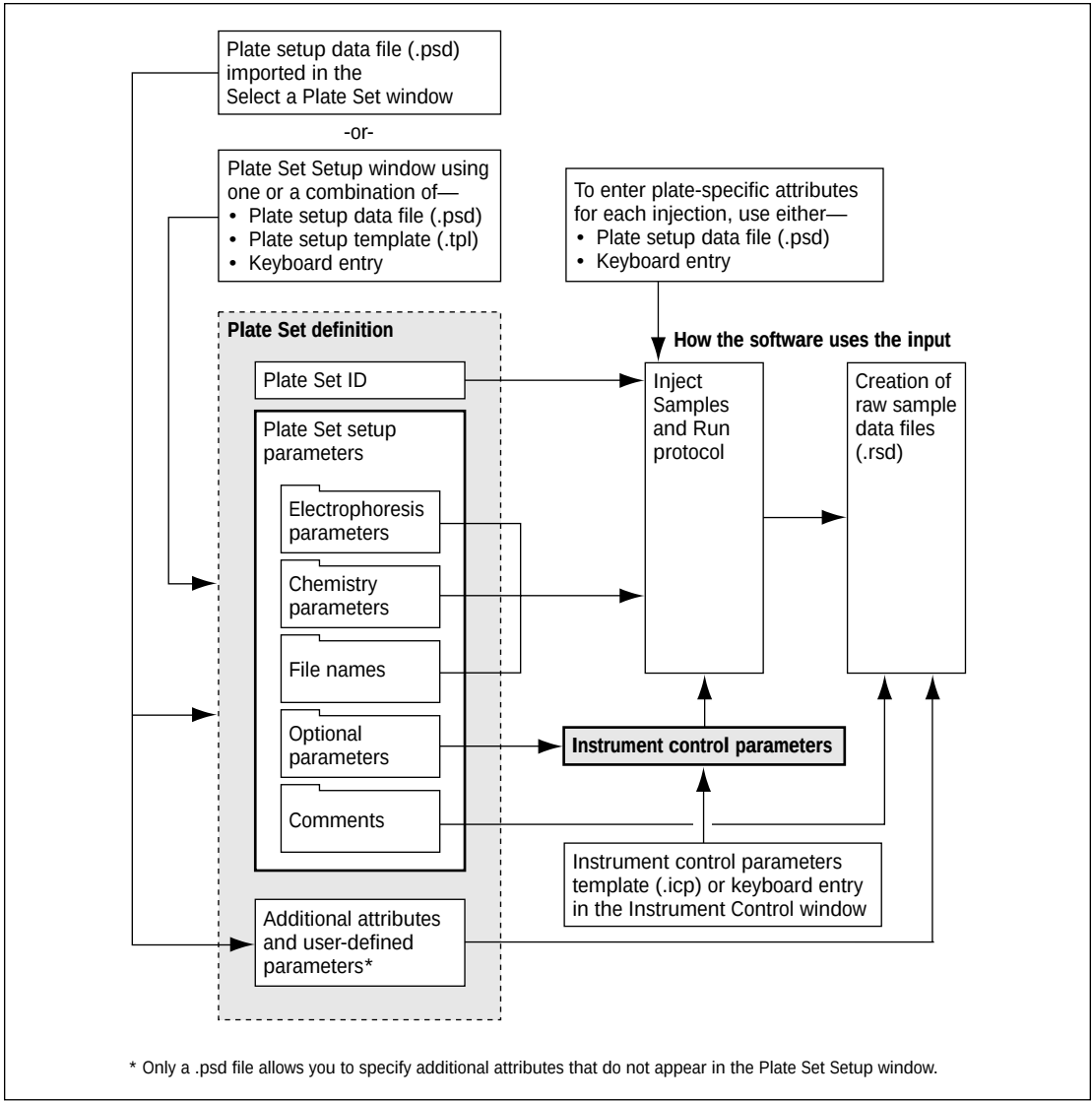


Figure 1-6. Overview of how the software uses the parameters specified for a SNP genotyping run.

1.6 File format overview

The initialization files (MegaBACE.ini, Chemistry.ini, and Basecall.ini), the parameter template files (.tpl and .icp), and the .psd files all share the same format. The file format consists of two columns (field column and value column) separated by one or more tabs. The software ignores additional character spaces. The order of the fields is arbitrary, and the field names are not case-sensitive. The Instrument Control Manager ignores any unrecognized fields or any blank rows between fields.

If you type a value that is invalid (for example an alpha character when the software requires a number, or a value that is out of range), the Instrument Control Manager displays an error message when the software loads the file. If this occurs, you must correct the value in the file before it can be loaded successfully.

1.7 File storage overview

The Instrument Control Manager looks for each type of file in a specific location—

- The MegaBACE.ini and Chemistry.ini files must be in the MegaBACE\Instrument Control\DataSystem folder.
- (Sequencing only) The Basecall.ini file must be in the MegaBACE\Base Calling folder.
- The template files (.tpl and .icp) must be stored in the MegaBACE\Instrument Control\Templates folder.
- The .psd files must be stored in the MegaBACE\Instrument Control\Psd folder or a location that you specify in the MegaBACE.ini file.

See chapter 15 for details on file storage.

Important

The system initialization files must be stored in the appropriate MegaBACE folder. The software will not function properly unless these files are stored properly. See chapter 15 for details.

Chapter 2 Specifying basic system preferences

This chapter describes the system initialization file that you can modify to configure the software for your laboratory. The topics are—

- Using the MegaBACE.ini file (section 2.1)
- Specifying a default application (section 2.2)
- Enabling or disabling the edit mode (section 2.3)
- Specifying how plate setup data (.psd) files are used (section 2.4)
- Using automated bar coding (section 2.5)
- Suppressing file creation for empty or bad capillaries (section 2.6)

2.1 Using the MegaBACE.ini file

The MegaBACE.ini file determines the initial settings the Instrument Control Manager uses to control the instrument. The Instrument Control Manager uses only one MegaBACE.ini file that loads automatically when you start the software. Table 2-1 describes the sections in the MegaBACE.ini file, and figure 2-1 shows the MegaBACE.ini file.

Table 2-1. Sections in the MegaBACE.ini file

Section name	Description	Reference
[General]	Lists the basic settings the Instrument Control Manager uses to run the instrument. The fields are active for all applications (sequencing, STR genotyping, and SNP genotyping).	Sections 2.1 through 2.5
[Sequencing]	Lists the default plate setup and instrument control templates for sequencing and all the base calling and file export modes.	Chapter 5
[Genotyping]	Lists the default plate setup and instrument control templates for STR genotyping.	Chapter 8
[SNP]	Lists the default plate set setup and instrument control templates for SNP genotyping and preferences for the injections during a multi-injection run.	Chapter 11
[Bad Capillaries]	Lists the well IDs of any bad capillaries (broken or clogged). The Instrument Control Manager does not create data files for capillaries that are empty, broken, or clogged. Note: You can specify empty capillaries in a plate setup template (.tpl) or a plate setup data file (.psd).	Section 2.6

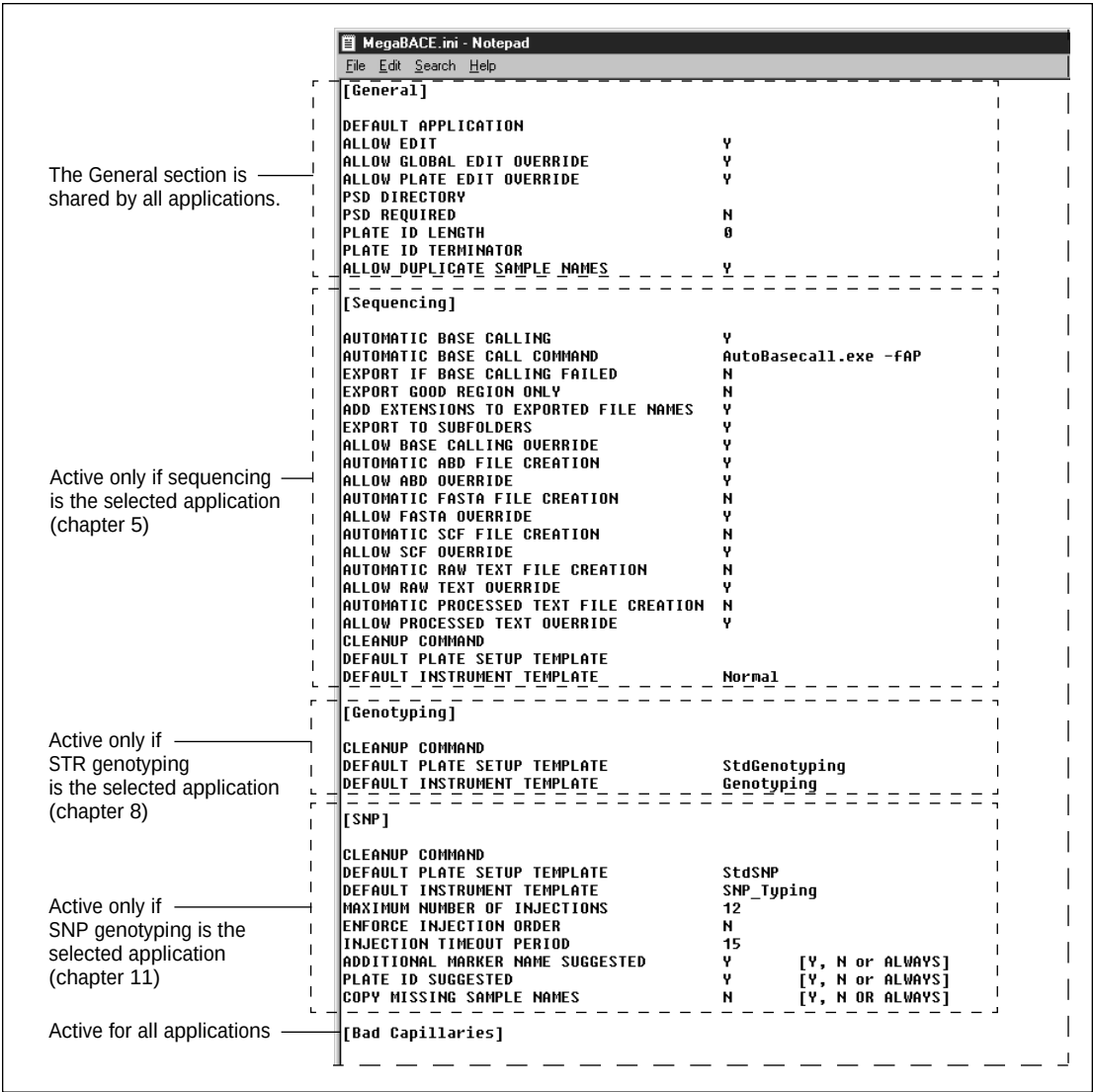


Figure 2-1. The MegaBACE.ini file opened in Notepad.

Important

2.1.1 Modifying the MegaBACE.ini file

The Instrument Control Manager can use the MegaBACE.ini file only if the file name is MegaBACE.ini and the file is stored in the ...\\MegaBACE\\Instrument Control\\DataSystem folder.

The MegaBACE.ini file is in a two-column, tab-separated format (figure 2-1). The first column (field column) contains the various modes and conditions that you can customize for the laboratory. The second column (value column) defines the setting. See section 1.6 for additional details on file formats.

To modify the MegaBACE.ini file—

1. Open the **MegaBACE.ini file** in the ...\\MegaBACE\\Instrument Control\\DataSystem folder using any text editor. See figure 2-1 for an example of the MegaBACE.ini file opened in Notepad.
2. Change the values in the value column as needed. Each field in the file has a specified data type. The data types for the value column include—
 - **Y/N**—You type **Y** for yes or **N** for no to set the conditions.
 - **String**—You type text.
 - **Integer**—You type a number.

Table 2-2 provides a description of the field names and the associated data type required for each field in the General section of the MegaBACE.ini file. For details on the application-specific sections of the MegaBACE.ini file, see the corresponding chapter for each application (chapter 5 for sequencing, chapter 8 for STR genotyping, or chapter 11 for SNP genotyping).

3. Save the file in the ...\\MegaBACE\\Instrument Control\\DataSystem folder.

Important

You must restart the Instrument Control Manager in order for the changes you make in the MegaBACE.ini file to take effect.

2.1.2 General field names

Table 2-2. The general field names and default settings of the MegaBACE.ini file

Field name	Description	Data type	Default
DEFAULT APPLICATION	Specifies the application to use as the default when an operator opens the Instrument Control Manager.	String [†]	Sequencing is the default.
ALLOW EDIT	Determines whether an operator can edit the parameters in the Plate Setup window and Instrument Control window. Y = operator can change parameters (edit mode is turned on) N = operator cannot change parameters (edit mode is turned off)	Y or N	Y
ALLOW GLOBAL EDIT OVERRIDE	Determines whether an operator can change the Edit Mode command on the Configure menu at any time. Y = operator can change mode N = operator cannot change mode	Y or N	Y
ALLOW PLATE EDIT OVERRIDE	Determines whether an operator can change the Edit Mode command on the Configure menu only for the plate definition that is currently being created (if Allow Global Edit Override is N). Y = operator can change mode N = operator cannot change mode	Y or N	Y
PSD DIRECTORY	Specifies the path for the psd folder.	Path [†]	...\MegaBACE\Instrument Control\Psd
PSD REQUIRED	Determines whether an operator can proceed without a .psd file (see section 2.4.2). Y = operator cannot proceed without a .psd file N = to create a plate definition, the operator can choose whether to import a .psd file	Y or N	Y
PLATE ID LENGTH	Determines the number of characters allowed in the bar code and eliminates the need to press ENTER.	Integer [‡]	N/A

Table 2-2. The general field names and default settings of the MegaBACE.ini file (continued)

Field name	Description	Data type	Default
PLATE ID TERMINATOR	Determines what terminates the entry of the plate ID, for example a character or group of characters. Eliminates the need to press ENTER (section 2.5.2).	String	None
ALLOW DUPLICATE SAMPLE NAMES	Allows an operator to use duplicate sample names. Y = allow duplicate names (Well IDs are appended to the sample names to create unique file names.) N = do not allow duplicate names	Y or N	Y

* String = A text string
† Path = The complete path to a file folder
‡ Integer = A number

2.2 Specifying a default application

Whenever you restart the Instrument Control Manager, the software automatically resets the selected application (sequencing, STR genotyping, or SNP genotyping) to the default specified in the MegaBACE.ini file. After the software starts, the operator can change the application selected in the Instrument Control Manager by choosing the Applications command from the Configure menu. The name of the selected application appears in the title bar of each of the Instrument Control Manager windows.

The Instrument Control Manager windows display only the parameters that are applicable for the selected application. For example, if sequencing is selected, the Plate Setup window displays only the sequencing chemistry names, and the Plate Catalog displays only the sequencing plate IDs.

To set the default application, you use the following field in the MegaBACE.ini file:

Field in MegaBACE.ini file	Value	Results
DEFAULT APPLICATION	(blank)	The sequencing application is the selected default.
	Genotyping	The STR genotyping application is the selected default.
	SNP Genotyping	The SNP genotyping application is the selected default.

For each application, you can specify the default plate setup template and the default instrument control template. For details on the application-specific sections of the MegaBACE.ini file, see the chapter for each application (chapter 5 for sequencing, chapter 8 for STR genotyping, or chapter 11 for SNP genotyping).

2.3 Enabling or disabling the edit mode

The edit mode determines whether an operator can manually edit the parameters displayed in the Plate Setup window or Plate Set Setup window and the Instrument Control window. The edit mode also determines whether an operator can change templates or select templates. If the edit mode is turned on, a check mark appears next to the Edit Mode command on the Configure menu. If the edit mode is turned off, no check mark appears.

The MegaBACE.ini file contains the default settings for the edit mode conditions (table 2-3). You can use the MegaBACE.ini file to enable or disable the edit mode menu command, which allows an operator to turn on or off the edit mode using the Configure menu. You can also specify whether the edit mode remains on or applies to one plate only.

Table 2-3 provides an overview of how to specify the edit mode conditions in the MegaBACE.ini file and how the Edit Mode command appears in the Configure menu.

Table 2-3. The edit mode fields in the MegaBACE.ini file

Allow Edit	Allow Global Edit Override	Allow Plate Edit Override	Edit Mode command in the Configure menu
Y	N/A	N/A	Turned on (checked). Can be turned off.
N	Y	N/A	Turned off (unchecked). Can be turned on at will.
N	N	Y	Turned off (unchecked). Can be turned on for one newly created plate only, and turns off automatically for all other plates.
N	N	N	Turned off (unchecked). No editing is allowed. An operator can create a new plate definition only by importing a .psd file or by using the default plate setup template.

The MegaBACE.ini file contains the following fields that allow you to specify the edit mode conditions:

- **Allow Edit**—Specifying **Y** turns the edit mode on and allows an operator to edit any of the plate setup or instrument parameters. The Edit Mode command on the Configure menu has a check mark to indicate that it is turned on. Specifying **N** turns off the edit mode, but an operator can turn on the edit mode if a **Y** is specified for the ALLOW GLOBAL EDIT OVERRIDE field.
- **Allow Global Edit Override**—Specifying **Y** allows an operator to turn on or off the edit mode in the Configure menu at will. While in the edit mode, an operator can edit the plate setup parameters or instrument parameters. Specifying **N** turns off the edit mode, which prevents the operator from turning on the edit mode, unless a **Y** is specified for the ALLOW PLATE EDIT OVERRIDE field.
- **Allow Plate Edit Override**—Specifying **Y** allows an operator to choose the edit mode for the current plate only (newly created plate). The edit mode becomes unavailable when the operator selects another plate. Specifying **N** disables the edit mode feature.

If all the above conditions are set to **N** in the MegaBACE.ini file, an operator cannot edit any of the plate setup parameters (plate set setup for SNP genotyping) or instrument parameters in the Instrument Control Manager windows. In this case, an operator can set up plate definitions (plate set definitions for SNP genotyping) only by importing .psd files or by using the default plate setup template.

When you turn on the edit mode, an operator can create a plate setup template, or select or set a default template from the Plate Setup Template menu. An operator can also create an instrument parameter template, or select a template from the Instrument Templates menu.

2.4 Specifying how plate setup data (.psd) files are used

A plate setup data file (.psd) allows you to define—

- **Sequencing and STR genotyping**—A plate ID, the plate setup parameters, and other plate- and well-level attributes so that an operator can automatically import the plate definition. The operator can import the .psd file either before the run in the Plate Setup window or at the start of the run in the Select a Plate window (*MegaBACE 4000 Instrument User's Guide*).



- **SNP genotyping**—A plate set ID, the plate set setup parameters, and other plate- and well-level attributes so that an operator can automatically import the plate set definition and plate-attributes for each injection. The operator can import the master .psd file containing the plate set definition and the attributes for the first plate either before the run in the Plate Set Setup window or at the start of the run in the Select a Plate Set window (*MegaBACE 4000 Instrument User's Guide*). The operator can import a plate-only .psd file for each subsequent injection.

With automatic import, the plate setup parameters or plate set setup parameters and plate attributes are specified with less operator intervention. The parameters that you specify in the .psd file become part of the extended header of each raw sample data file (.rsd) created during the run.

2.4.1 Specifying a new .psd file folder location

By default, the .psd files are saved in the ...\\MegaBACE\\Instrument Control\\Psd folder. To specify a new storage location for the .psd files, you use the following field in the MegaBACE.ini file:

Field in MegaBACE.ini file	Value
PSD DIRECTORY	...\\MegaBACE\\Instrument Control\\Psd

You can specify a file path on a network so that the .psd files can be shared by multiple workstations. To do this, use the full network path location. For example—

Field in MegaBACE.ini file	Value
PSD DIRECTORY	\\computer_name\\folder_name\\Psd

2.4.2 Specifying that .psd files are required

The operator can always import a .psd file. However, the PSD REQUIRED field determines whether the operator can proceed without a .psd file. By default, the value for this field is N. If the PSD REQUIRED field is set to Y, the operators must use a .psd file if they use the Plate Setup window to create plate definitions. If the software cannot find or load the .psd file, an error message appears and the operation stops.



If the PSD REQUIRED field is set to N, the operator is not required to use a .psd file when creating the plate definitions in the Plate Setup window. The software still accepts .psd files.

Field in MegaBACE.ini file	Value
PSD REQUIRED	Y or N

2.5 Using automated bar coding

Important

Regardless of the method used to enter the plate ID, the final text entry or bar code (less a possible plate ID terminator) must match the name (case insensitive) of the .psd file in order for a .psd file to be imported.

To eliminate the need to press ENTER when importing a .psd file, you can use the following fields in the MegaBACE.ini file:

- PLATE ID LENGTH
- PLATE ID TERMINATOR

(Sequencing and STR genotyping) When an operator types a .psd file name or scans a bar code in the Plate ID box of the Plate Setup window or the Select a Plate window, the software automatically imports the .psd file when the PLATE ID LENGTH or PLATE ID TERMINATOR conditions are met.

(SNP genotyping) When an operator types a .psd file name or scans a bar code in the Plate Set Setup window, the Select a Plate Set window, or the Bar Code box of the Sample Injection window, the software automatically imports the .psd file when the PLATE ID LENGTH or PLATE ID TERMINATOR conditions are met.

2.5.1 Specifying the plate ID length

You can predefine the number of characters in the bar code or .psd file name using the PLATE ID LENGTH field.

Field in MegaBACE.ini file	Value
PLATE ID LENGTH	Integer

When the number of characters entered reaches the predefined value, the software automatically imports the .psd file.



2.5.2 Specifying the plate ID terminator

You can determine which character or group of characters terminates the entry of the bar code or .psd file name. You can use any printable character as the plate ID terminator. After detecting the terminator in the string entered, the Instrument Control Manager stops accepting characters, removes the terminator from the end of the string, and automatically imports the .psd file.

Field in MegaBACE.ini file	Value
PLATE ID TERMINATOR	String

2.6 Suppressing file creation for empty or bad capillaries

The Instrument Control Manager allows you to suppress the creation of raw sample data files (.rsd) for capillaries that are empty, broken, or clogged. You can list the well ID for each defective capillary in the [Bad Capillaries] section of the MegaBACE.ini file. This is a global setting. For every capillary listed, a red arrow appears above the capillary button in the Run Image window (figure 2-2).

You can specify empty capillaries in a .psd file, a .tpl file, or manually in the Plate Setup window (sequencing and STR genotyping) or the Plate Set Setup window (SNP genotyping). By specifying NOT USED as the file name, an empty capillary appears as black, and a yellow arrow appears above the well button in the Run Image window.



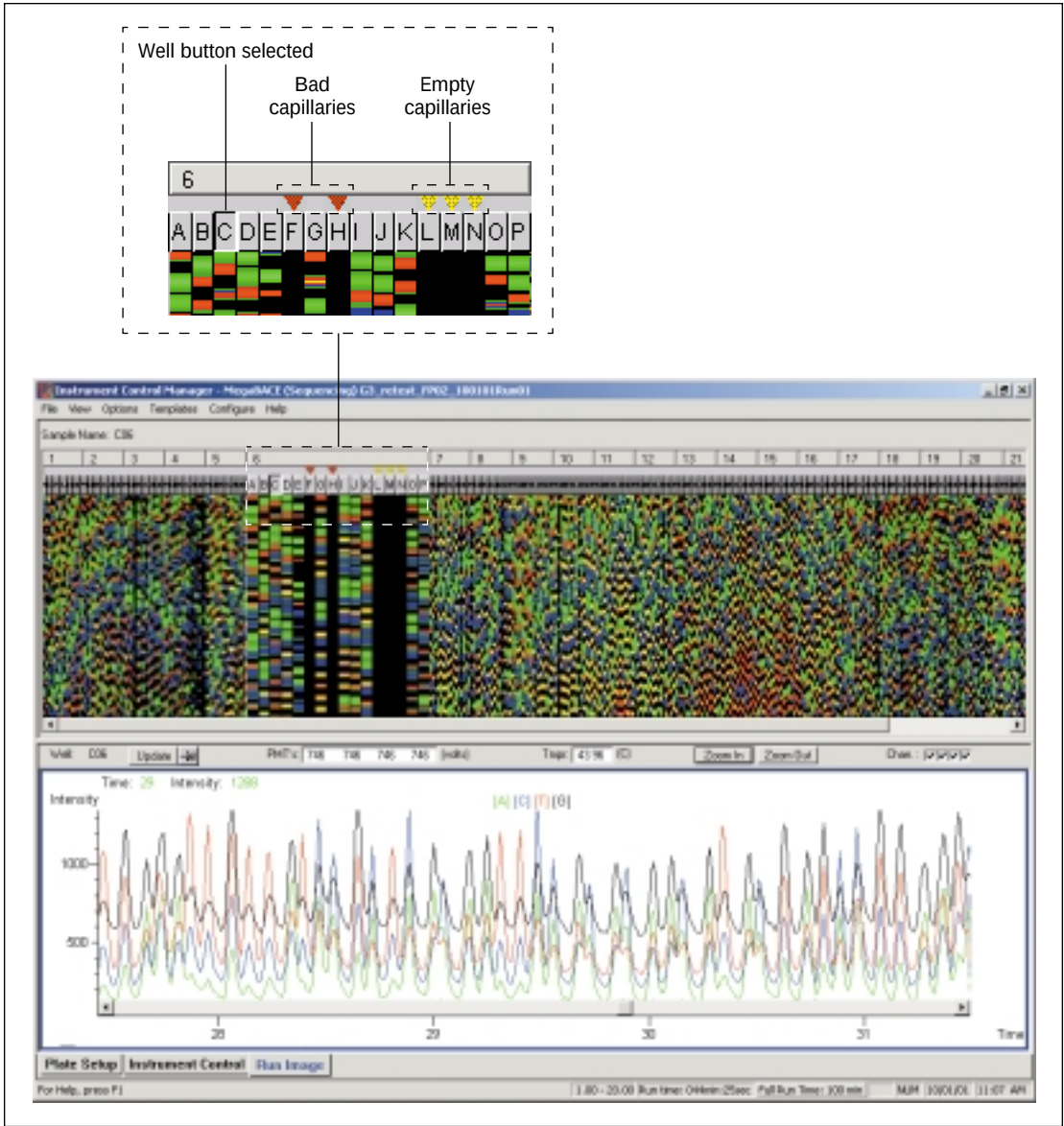


Figure 2-2. The Run Image window for an instrument with bad capillaries and empty capillaries indicated.

Chapter 3 Specifying the chemistry parameters

This chapter describes the system chemistry initialization file that you can modify to configure the system for your laboratory workflow. The topics are—

- How the software uses the Chemistry.ini file (section 3.1)
- About specifying the chemistry (section 3.2)
- The fields in the Chemistry.ini file (section 3.3)
- Adding a new set of chemistry parameters (section 3.4)

3.1 How the software uses the Chemistry.ini file

The Chemistry.ini file contains all the chemistry parameter sets that are available in the Instrument Control Manager. The Instrument Control Manager uses only one Chemistry.ini file to define all the available chemistry parameter sets. See section 1.6 for a description of the file format.

3.2 About specifying the chemistry

To use one of the available chemistry parameter sets, you must specify a chemistry name that is recognized by the Instrument Control Manager. You can specify the chemistry name in a—

- Plate setup data file (.psd)

Note: For SNP genotyping, if you are using a .psd file to specify the chemistry name, it must be a master .psd file.

- Plate setup template file (.tpl)

If you turn on the edit mode, an operator can select the chemistry name in the Plate Setup window (sequencing and STR genotyping) or the Plate Set Setup window (SNP genotyping). Only the names of the chemistry parameter sets for the selected application appear in the Chemistry Name list in the Chemistry Parameters tab (figure 3-1).

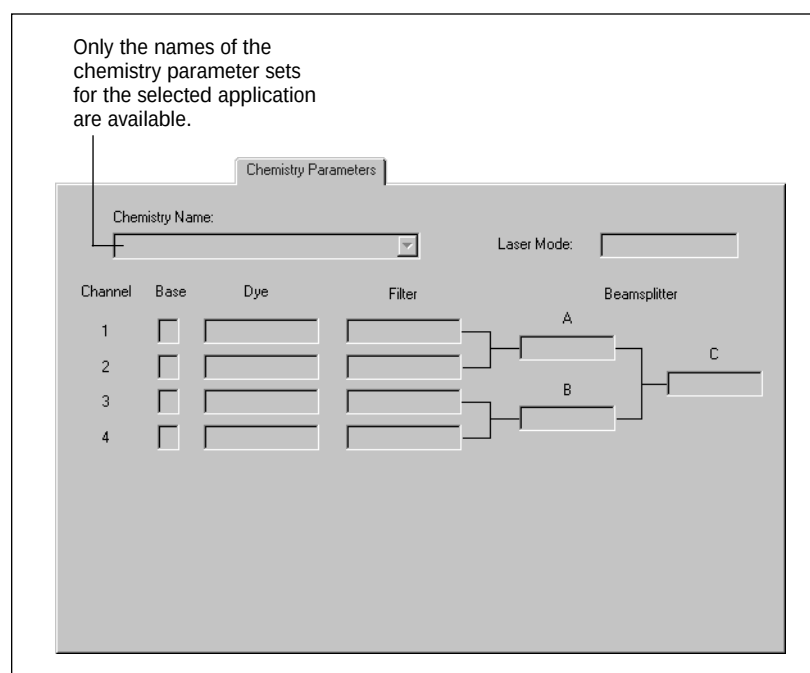


Figure 3-1. The Chemistry Parameters tab of the Plate Setup window (sequencing and STR genotyping) or the Plate Set Setup window (SNP genotyping).

3.3 The fields in the Chemistry.ini file

The fields in the Chemistry.ini file include the following:

- **Application**—Required field if you are defining a parameter set for an application other than sequencing. If you do not specify the application, the Instrument Control Manager assumes that the chemistry parameter set is for the sequencing application. The Instrument Control Manager recognizes the following three application types:
 - Sequencing
 - Genotyping (STR genotyping)
 - SNP genotyping
- **Beamsplitters**—BEAMSPLITTER A indicates the beamsplitter used in spectral channels 1 and 2. BEAMSPLITTER B indicates the beamsplitter used in spectral channels 3 and 4. BEAMSPLITTER C splits the emission light between BEAMSPLITTER A and BEAMSPLITTER B.

- **Filters**—The filter number indicates the spectral channel number. For example, FILTER 1 indicates the filter used to detect the emissions in spectral channel 1.
- **Bases**—Required field for the sequencing and SNP genotyping applications. The base number indicates the channel number used to detect the base. For example, BASE 1 indicates the base detected through spectral channel 1. Valid values are A, G, C, or T.
- **Dyes**—Required field for the STR genotyping application. The dye number indicates the channel number used to detect the dye. For example, DYE 1 indicates the dye detected through spectral channel 1.
- **Laser mode**—Required field. The only valid value for the MegaBACE 4000 instrument is blue.

3.4 Adding a new set of chemistry parameters

Important

If the application is not specified in the chemistry parameter set, the Instrument Control Manager assumes the chemistry is for the sequencing application, and you must specify the bases.

To add an additional chemistry parameter set to the Chemistry.ini file—

1. Open the Chemistry.ini file using a text editor or Excel. The Chemistry.ini file is located in the ...\\MegaBACE\\Instrument Control\\DataSystem folder. Figure 3-2 shows a section of the Chemistry.ini file.

The name of each chemistry parameter set appears in brackets, for example [TSII-Version 2 Terminators]. The field names appear in the left column, and the values appear in the right column. The columns are separated by one or more tabs.

2. Scroll to the end of the file and type a name for the new chemistry parameter set. You must provide a unique name for each chemistry parameter set. Make sure you put the name in brackets (figure 3-2).
3. Copy and paste the field names and values from a chemistry parameter set that is similar to the one you are adding. Make any required changes to the field values for the new parameter set.
4. Save the file. The Instrument Control Manager can use the file only if the file name is Chemistry.ini and the file is stored in the ...\\MegaBACE\\Instrument Control\\DataSystem folder.

Note: The changes appear immediately in the Instrument Control Manager without restarting the software.

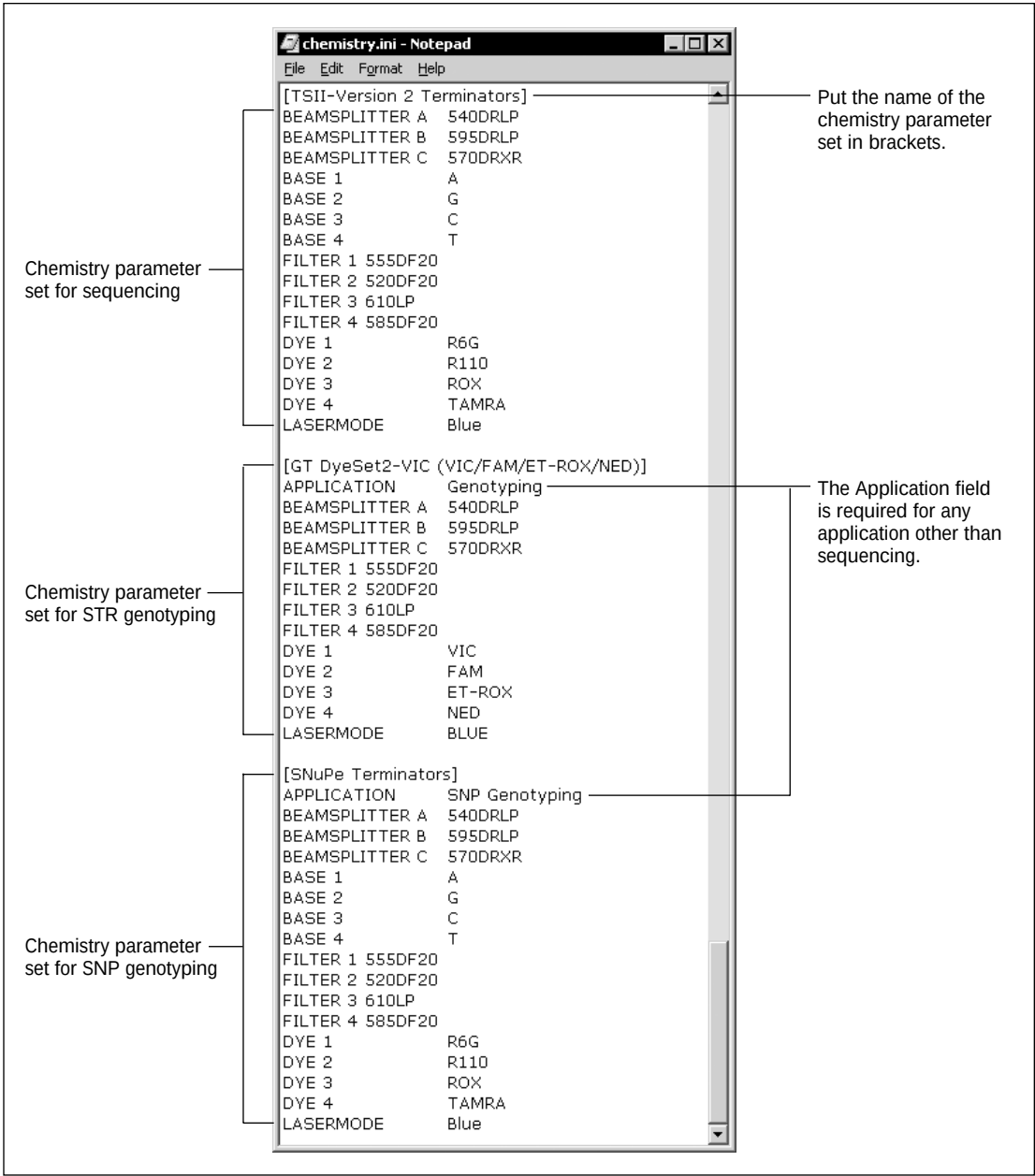


Figure 3-2. A section of the Chemistry.ini file opened in Notepad.

- 
5. Verify that the changes appear in the Chemistry Parameters tab (figure 3-1). If the name for a new or revised chemistry parameter set does not appear in the Chemistry Name list for the selected application, check the parameters in the Chemistry.ini file for syntax errors.

Syntax errors can include—

- The required field is missing. The required fields for sequencing are laser mode and base names.
- The required field has an out-of-range value. For example, a base name has a letter other than A, G, C, or T for the value.
- Duplicate values exist for different fields, for example, two bases are named A.
- The parameter set includes a misspelled field name or value for a required field.

Chapter 4 Creating templates for the instrument control parameters

The Instrument Control Manager allows you to save the frequently used instrument control parameters as templates so that you do not have to reenter the parameters. This chapter describes how to create the instrument control templates. The topics are—

- About instrument control templates (section 4.1)
- Using the Instrument Control window to create a template (section 4.2)
- Using a text file to create an instrument control template (.icp) (section 4.3)

4.1 About instrument control templates

The instrument control parameters are displayed in the Run Parameters tab and Sleep Parameters tab of the Instrument Control window (figure 4-1). The software includes the following instrument control template files (.icp):

- **Sequencing**—The Normal.icp template is provided.
- **STR genotyping**—The Genotyping.icp template is provided.
- **SNP genotyping**—The SNP_Typing.icp template is provided.

You can use the templates provided, or modify the existing templates to create new templates. You can also create a template from the values entered in the Instrument Control window, or you can create a template using a text editor or Excel. Tables 4-1 and 4-2 describe the instrument control parameters.

You can specify a default .icp file in the MegaBACE.ini file for each application (chapter 5 for sequencing, chapter 8 for STR genotyping, or chapter 11 for SNP genotyping). The default .icp file is loaded into the Instrument Control window at startup. If you enable the edit mode (section 2.3), an operator can also select a template from the Instrument Templates menu.

Part one General system configuration

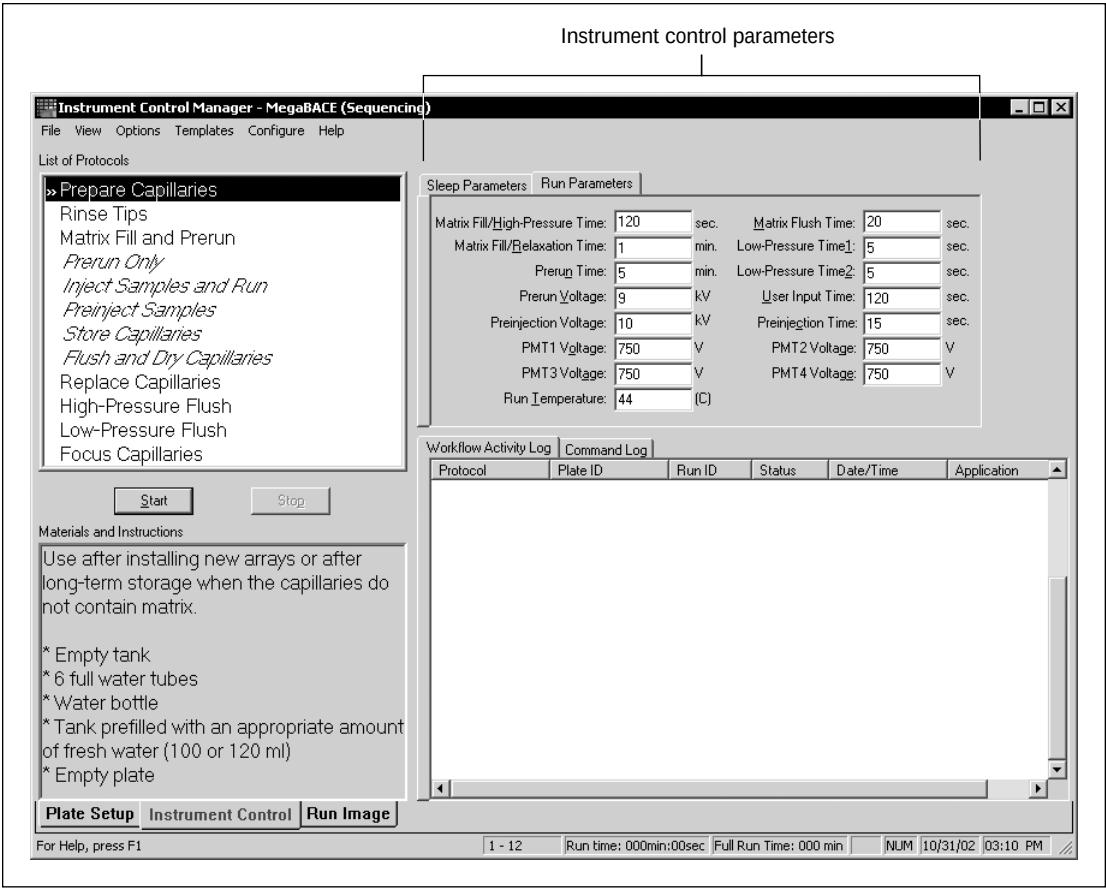


Figure 4-1. The Instrument Control window with the sequencing instrument parameters displayed.

Table 4-1. Instrument control run parameters

Field name	Description (value range)
Matrix Fill/High-Pressure Time	The time (seconds) for applying high pressure to fill the capillaries with matrix (range: 1–600 s).
Matrix Fill/Relaxation Time	The time (minutes) to allow the matrix to equilibrate (range: 1–120 min).
Prerun Time	The time (minutes) for the electrophoresis prerun (range: 1–120 min).
Prerun Voltage	The voltage for the prerun (range: 1–15 kV).
Preinjection Voltage	The voltage for preinjecting samples (range: 1–15 kV).
PMT1 Voltage	The voltage for PMT1 (range: 500–900 V).
PMT2 Voltage	The voltage for PMT2 (range: 500–900 V).
PMT3 Voltage	The voltage for PMT3 (range: 500–900 V).
PMT4 Voltage	The voltage for PMT4 (range: 500–900 V).
Run Temperature	The electrophoresis compartment temperature for the run (range 27–44 °C).
Matrix Flush Time	The time (seconds) for the first matrix flush (range: 1–60 s).
Low-Pressure Time1	The time (seconds) for the nitrogen pressure to flush water through the capillaries (range: 1–60 s).
Low-Pressure Time2	The time (seconds) for the nitrogen pressure to dry the capillaries (range: 1–300 s).
User Input Time	The time (seconds) you need to open and close the anode and cathode drawers after the display on the instrument gives the instruction to load the plate, tank, or tubes in either the anode or the cathode side of the instrument (range: 10–240 s).
Preinjection Time	The time (seconds) for preinjecting samples (range: 1–600 s).

Table 4-2. Instrument control sleep parameters

Field name	Description and value range
Sleep After This Run	This setting is not saved in the template. The check box that you check if you want the instrument set in sleep mode after the current or next run.
Sleep Temperature	Although the Sleep Temperature field requires a value, the software turns off the electrophoresis compartment heater during the Sleep After This Run and the Store Capillaries protocols. The value range is 25–30.
Sleep Time	The length of time (hours) the instrument will be left idle, during which the software turns off the electrophoresis compartment heater and the laser. The range is 1–168 h (Store Capillaries protocol) or 1–16 h (Sleep After This Run option).

4.2 Using the Instrument Control window to create a template

The MegaBACE software includes instrument control templates (.icp files) that you can use. You can modify the existing templates to create new templates. You can create a template from the values entered in the Instrument Control window (figure 4-1), or you can create a template using a text editor or Excel (section 4.3). See tables 4-1 and 4-2 for a list of the parameters.

Important The edit mode must be turned on to allow creating, editing, saving, or importing instrument control templates in the Instrument Control Manager. See section 2.3 for information about enabling the edit mode in the MegaBACE.ini file.

To create an instrument control template—

1. Start the Instrument Control Manager. In the **MegaBACE** folder on the Windows desktop, double-click the **Instrument Control Manager** icon.
2. Click the **Instrument Control** tab to display the Instrument Control window (figure 4-1).
3. To modify an existing template—
 - On the **Templates** menu, point to **Instrument Template**, and click **Select Template**. The Select Template window appears and displays the list of available templates (.icp).

- Select the template (.icp) you want to modify, and click **Open**. The values of the selected template appear in the Instrument Parameters area (figure 4-1).
4. Modify any of the parameters.
 5. To save the current instrument control parameters as a template—
 - On the **Templates** menu, point to **Instrument Template**, and click **Save Template**. The Save window appears.
 - Type a new name, and click **Save**. The template is saved in the ...\\MegaBACE\\Instrument Control\\Templates folder. The software adds the .icp file extension.

4.3 Using a text file to create an instrument control template (.icp)

You can also create an .icp file in a two-column format in Excel or in any text editor using tabs as column separators. You must use the field names that the Instrument Control Manager recognizes and type them in the first column (field column) as shown in figure 4-2. You type the values for each field in the second column (value column). The Instrument Control Manager uses the values only for the fields it recognizes. See table 4-1 for a list of the run parameter (field) names and a description of the expected values. See table 4-2 for a list of sleep parameters.

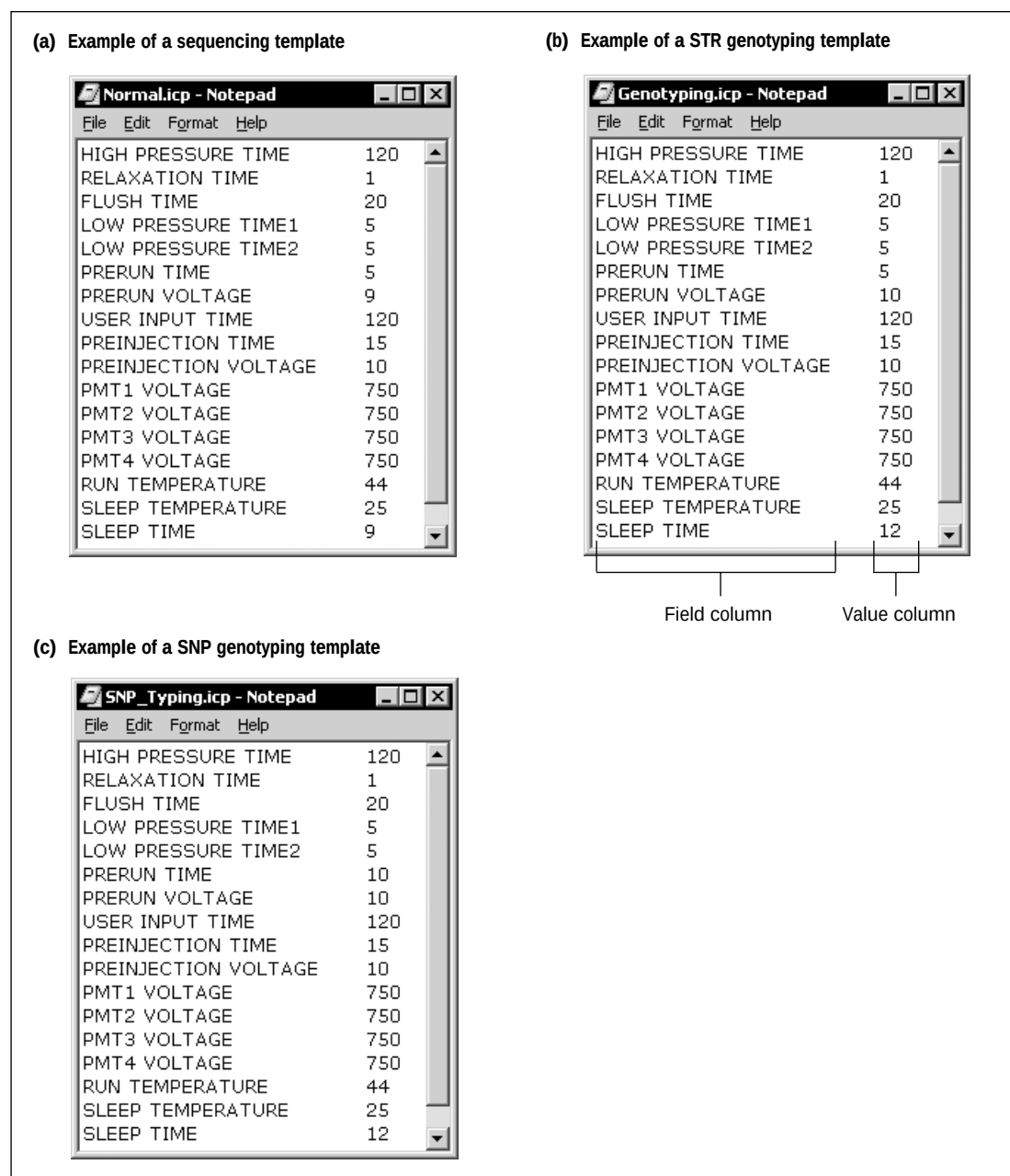


Figure 4-2. Examples of instrument control template files (.icp) created using Notepad.

Part two

Setting up the system for sequencing

Chapter 5 Specifying preferences for sequencing runs

This chapter describes how to set up the automatic base calling feature and how to automatically export files to other file formats. The topics are—

- The sequencing fields in the MegaBACE.ini file (section 5.1)
- Setting default parameter templates (section 5.2)
- Configuring the automatic base calling feature (section 5.3)
- Turning on or off the automatic base calling feature (section 5.4)
- Setting up the automatic file export feature (section 5.5)
- Settings in the Basecall.ini file (section 5.6)
- Specifying which base caller to use (section 5.7)

Make sure you configure the General section of the MegaBACE.ini file (chapter 2) for your laboratory. To create additional chemistry parameters, see chapter 3.

5.1 The sequencing fields in the MegaBACE.ini file

The MegaBACE.ini file determines the settings available in the Instrument Control Manager when you open the software. The file is divided into sections: General, Sequencing, Genotyping (STR genotyping), and SNP (SNP genotyping). Figure 5-1 shows the MegaBACE.ini file with the Sequencing fields highlighted.

Table 5-1 lists the field names and default settings in the Sequencing section of the MegaBACE.ini file. See chapter 2 for a description of the other fields in the MegaBACE.ini file.

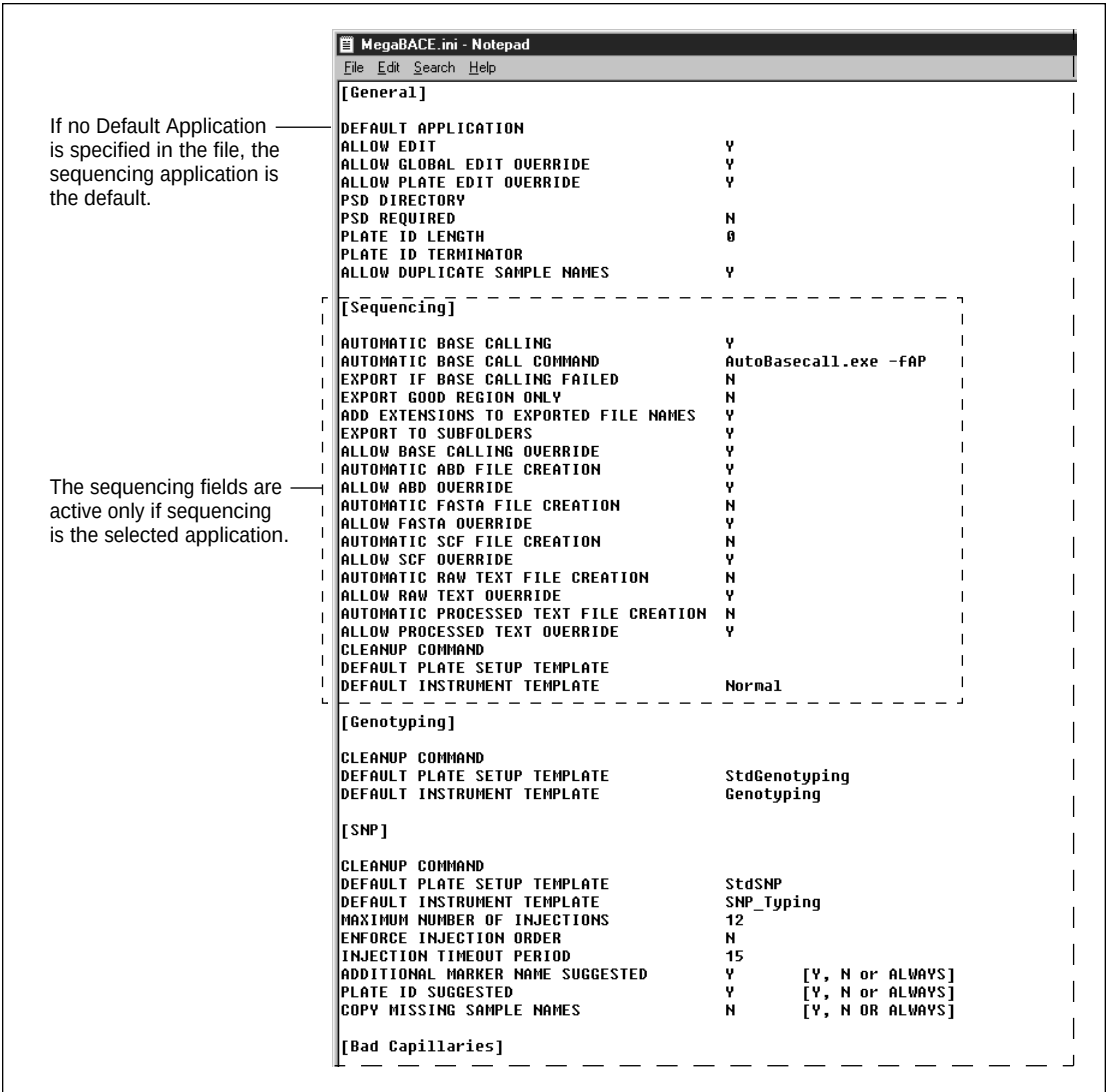


Figure 5-1. The MegaBACE.ini file with the sequencing fields highlighted.

Table 5-1. The field names and default settings for the sequencing section of the MegaBACE.ini file

Field name	Description	Setting type	Default
AUTOMATIC BASE CALLING	Determines whether the Instrument Control Manager launches automatic base calling. Y = launch base calling N = do not launch base calling	Y or N	Y
AUTOMATIC BASE CALL COMMAND	Allows automatic base calling to start right after data collection. Caution: If this field is blank, automatic base calling will fail. The field also allows you to specify a FASTA option: -fa, -fap, and -fag (section 5.5.2).	String	AutoBasecall.exe -fap
EXPORT IF BASE CALLING FAILED	Allows an operator to specify whether to export files for which base calling failed. Y = export files N = do not export files	Y or N	N
EXPORT GOOD REGION ONLY	Allows an operator to specify whether to export only the good region of sequence. Y = export good region N = export all	Y or N	N
ADD EXTENSIONS TO EXPORTED FILE NAMES	Adds standard extensions to exported file names. Y = add extensions N = do not add extensions	Y or N	Y
EXPORT TO SUBFOLDERS	Exports specified files to subfolders named with the file extensions. Y = use subfolders N = do not use subfolders	Y or N	Y
ALLOW BASE CALLING OVERRIDE	Determines whether an operator can override (turn on or turn off) automatic base calling from the Automatic Base Calling window, which opens from the Configure menu. Y = operator can override N = operator cannot override	Y or N	Y

Table 5-1. The field names and default settings for the sequencing section of the MegaBACE.ini file (continued)

Field name	Description	Setting type	Default
AUTOMATIC ABD FILE CREATION	Determines whether the Instrument Control Manager automatically creates ABD files after base calling. Y = create ABD files N = do not create ABD files	Y or N	Y
ALLOW ABD OVERRIDE	Determines whether an operator can override automatic ABD file creation from the Automatic Base Calling window, which opens from the Configure menu. Y = operator can override N = operator cannot override	Y or N	Y
AUTOMATIC FASTA FILE CREATION	Determines whether the Instrument Control Manager automatically creates FASTA files (.seq) after base calling. Y = create FASTA files N = do not create FASTA files	Y or N	N
ALLOW FASTA OVERRIDE	Determines whether an operator can override automatic FASTA file creation from the Automatic Base Calling window, which opens from the Configure menu. Y = operator can override N = operator cannot override	Y or N	Y
AUTOMATIC SCF FILE CREATION	Determines whether the Instrument Control Manager automatically creates SCF files after base calling. Y = create SCF files N = do not create SCF files	Y or N	N
ALLOW SCF OVERRIDE	Determines whether an operator can override automatic SCF file creation from the Automatic Base Calling window, which opens from the Configure menu. Y = operator can override N = operator cannot override	Y or N	Y

Table 5-1. The field names and default settings for the sequencing section of the MegaBACE.ini file (continued)

Field name	Description	Setting type	Default
AUTOMATIC RAW TEXT FILE CREATION	Determines whether the Instrument Control Manager automatically creates raw text files after base calling. Y = create raw text files N = do not create raw text files	Y or N	N
ALLOW RAW TEXT OVERRIDE	Determines whether an operator can override automatic raw text file creation from the Automatic Base Calling window, which opens from the Configure menu. Y = operator can override N = operator cannot override	Y or N	Y
AUTOMATIC PROCESSED TEXT FILE CREATION	Determines whether the Instrument Control Manager automatically creates processed text files after base calling. Y = create processed text files N = do not create processed text files	Y or N	N
ALLOW PROCESSED TEXT OVERRIDE	Determines whether an operator can override automatic processed text file creation from the Automatic Base Calling window, which opens from the Configure menu. Y = operator can override N = operator cannot override	Y or N	Y
CLEANUP COMMAND	Allows the administrator to specify a program or batch file to be executed at the end of the run after base calling. This can be useful for transferring plate records and raw data files to another workstation for analysis or archiving. See appendix C for a more detailed description.	String consisting of the path to the program or batch file (full or relative path to the instrument control folder)	None
DEFAULT PLATE SETUP TEMPLATE	Specifies a default template for the plate setup parameters.	String	None
DEFAULT INSTRUMENT TEMPLATE	Specifies a default template for the instrument control parameters.	String	Normal

5.2 Setting default parameter templates

In the MegaBACE.ini file, you can specify the—

- Default plate setup template
- Default instrument template

Whenever the Instrument Control Manager is restarted, the software automatically resets the templates selected for a given application to the default templates specified in the MegaBACE.ini file. However, if the edit mode is on, the operator can choose a different template in the Instrument Control Manager.

5.2.1 Setting the default plate setup template

If you specify a default plate setup template in the MegaBACE.ini file, the values from the plate setup template appear automatically in the Plate Setup window when the operator clicks **New**. To specify a default template for sequencing, use the following field in the sequencing section of the MegaBACE.ini file:

Field in MegaBACE.ini file	Value
DEFAULT PLATE SETUP TEMPLATE	Type the file name of the .tpl file, without the .tpl file extension.

Note: If a plate setup data file (.psd) is used to define the plate for a run, the values for any duplicate parameters overwrite the values specified on the screen which typically come from the default plate setup template. The original default values are automatically restored at the completion of the run.

5.2.2 Setting the default instrument control template

If you specify a default instrument control template in the MegaBACE.ini file, the values appear automatically in the Instrument Control window when the Instrument Control Manager starts. To specify a default template for sequencing, use the following field in the sequencing section of the MegaBACE.ini file:

Field in MegaBACE.ini file	Value
DEFAULT INSTRUMENT TEMPLATE	Type the file name of the .icp file, without the .icp file extension.

5.3 Configuring the automatic base calling feature

The software can perform base calling automatically after the data collection is completed for a run. To configure the automatic base calling feature, you use a combination of the following files:

- **MegaBACE.ini file**—The Sequencing section of this file allows you to specify if the automatic base calling feature is on and whether an operator can turn on or off the feature in the Instrument Control Manager (section 5.4). This file also allows you to configure the file export feature (section 5.5).
- **Basecall.ini file**—The Base Callers section of the file lists the available base callers, and the Default section of the file can specify a default base caller (sections 5.6 and 5.7.1).
- **Plate setup template (.tpl)**—If used, a template can specify the base caller to be used for multiple plates (section 5.7.2 and chapter 6).
- **Plate setup data file (.psd)**—If used, a .psd file can specify both plate-level and well-level base callers for a given plate (section 5.7.2 and chapter 7).

5.4 Turning on or off the automatic base calling feature

The Sequencing section of the MegaBACE.ini file allows you to turn on the automatic base calling mode or to allow an operator to turn on or off the automatic base calling from the Configure menu in the Instrument Control Manager. The settings in the MegaBACE.ini file are tied directly to the settings in the Automatic Base Calling window in the Instrument Control Manager. Table 5-2 provides an overview of how you specify the automatic base calling conditions in the MegaBACE.ini file and how the Automatic Base Calling command appears in the Configure menu and in the Automatic Base Calling window (figure 5-2). See table 5-1 for a full list of the fields in the Sequencing section of the MegaBACE.ini file.

Table 5-2. The automatic base calling mode conditions in the MegaBACE.ini file

AUTOMATIC BASE CALLING field	ALLOW BASE CALLING OVERRIDE field	Automatic Base Call command in the Configure menu and the Automatic Base Calling window
Y	Y	Turned on (checked). Can be turned off in the Automatic Base Calling window.
Y	N	Turned on (checked). Cannot be turned off.
N	Y	Turned off (unchecked). Can be turned on at will in the Automatic Base Calling window.
N	N	Turned off and unchecked. Cannot be turned on.

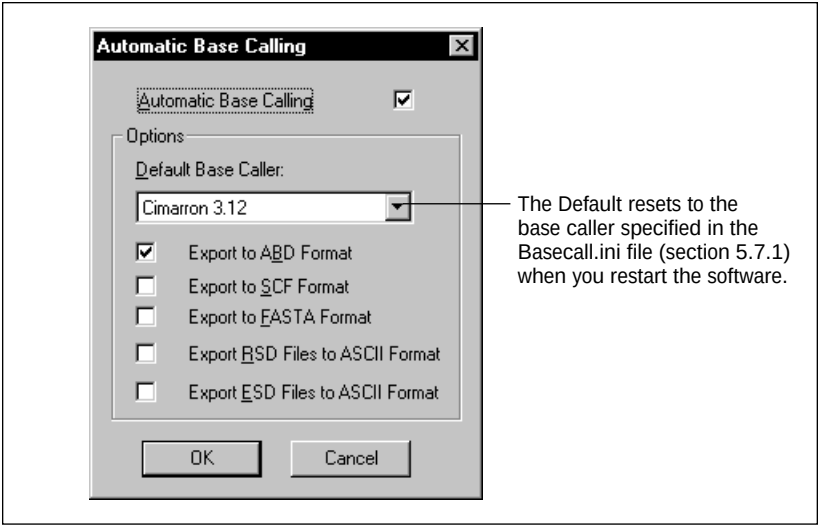


Figure 5-2. The Automatic Base Calling window, which is accessed from the Configure menu in the Instrument Control Manager.

Note: For a description of how to use the Automatic Base Calling window, see the *MegaBACE 4000 Instrument User's Guide*.

5.5 Setting up the automatic file export feature

After automatic base calling, the software can automatically export the resulting files to other file formats in addition to the base-called sample data files (.esd). You can also specify file export options, such as whether to export if base calling fails or whether to add extensions to the file names.

5.5.1 Specifying the file formats for export

The additional format options for automatic file export are—

- **ABD**—Contains all the raw data, including the electropherogram and spectral separation table. If base-called, the file also contains all the raw and base-called data, including the called sequence.
- **FASTA**—Contains the file name, file location, and the called sequence for each well in the base-called plate (default). See section 5.5.2.
- **SCF**—Contains the raw data and the electropherograms. If base-called, includes the raw and base-called data, and the called sequence. See section 5.5.3.
- **ASCII text**—Contains the plate name, sample name, instrument ID, base caller name, electrical current information, scan rate, instrument parameters. If base-called, includes the sequence starting and ending points and the number of bases called.

You can use the fields in the Sequencing section of the MegaBACE.ini file (table 5-1) to select the file formats you want the Instrument Control Manager to create. You also use the fields in the MegaBACE.ini file to specify whether the operator can override the MegaBACE.ini file export selections that appear in the Automatic Base Calling window (figure 5-2). In the MegaBACE.ini file, if you specify—

- **Y** for any of the export file formats, the corresponding check boxes for the specified file formats are selected in the Automatic Base Calling window.
- **N** for any of the export file formats, the corresponding check boxes for the specified file formats are not selected in the Automatic Base Calling window.
- **Y** to ALLOW [FILE_FORMAT] OVERRIDE, an operator can select or deselect the corresponding export options in the Automatic Base Calling window.

Note: Automatic base calling must be turned on to export the selected file format(s).

5.5.2 The FASTA file format

The FASTA file creation has three different options, and you can specify one, two, or all three options. If you want to use the FASTA export feature, you must select at least one of the following options:

- **-fa** for creation of one .seq file for each analyzed .esd file. These files are stored in the ...\\MegaBACE\\AnalyzedData\\plateID_runID folder.
- **-fap** for creation of one .seq file for each base-called plate. This file is stored in the ...\\MegaBACE\\AnalyzedData\\plateID_runID folder. This is the default switch used for FASTA in the Automatic Base Call Command in the MegaBACE.ini file.
- **-fag** for global creation of one .seq file for all the analyzed plates. This file is stored in the ...\\MegaBACE\\AnalyzedData\\plateID_runID folder.

To specify the FASTA file format, you use the AUTOMATIC BASE CALL COMMAND field in the MegaBACE.ini file. For example, to specify the -fap format, you would use the following:

Field in MegaBACE.ini file	Value
AUTOMATIC BASE CALL COMMAND	AutoBasecall.exe -fap

5.5.3 Changing the SCF version

For the SCF format, you can choose whether to export in SCF 2.00 or SCF 3.00 (default) file format. To change the SCF version, you use the Basecall.ini file (section 5.6).

5.5.4 Using subfolders for the exported files

You can use the EXPORT TO SUBFOLDERS field in the MegaBACE.ini file (table 5-1) to specify whether to export the selected file formats to subfolders. If you specify—

- **Y**—The software creates subfolders, within the raw or analyzed run folder, respectively, for the specified export files. Each subfolder has the name of the export file format (ABD, SCF, FASTA, ASCII text). See table 15-3 for a detailed description of each file format.
- **N**—The software exports the files in the selected formats to the top-level raw or analyzed run folder, respectively.

5.5.5 Using the file export options

If you turn on the automatic base calling feature, the Instrument Control Manager provides the following export options:

- **EXPORT IF BASE CALLING FAILED**—Default is **N**.
- **EXPORT GOOD REGION ONLY**—Default is **N**.
- **ADD EXTENSIONS TO EXPORTED FILE NAMES**—Default is **Y**.

You turn on or off the options in the MegaBACE.ini file by typing **Y** or **N**.

5.6 Settings in the Basecall.ini file

Figure 5-3 shows the Basecall.ini file. The Basecall.ini file is in a two-column, tab-separated format. Table 5-3 provides a description of each section and field in the Basecall.ini file.

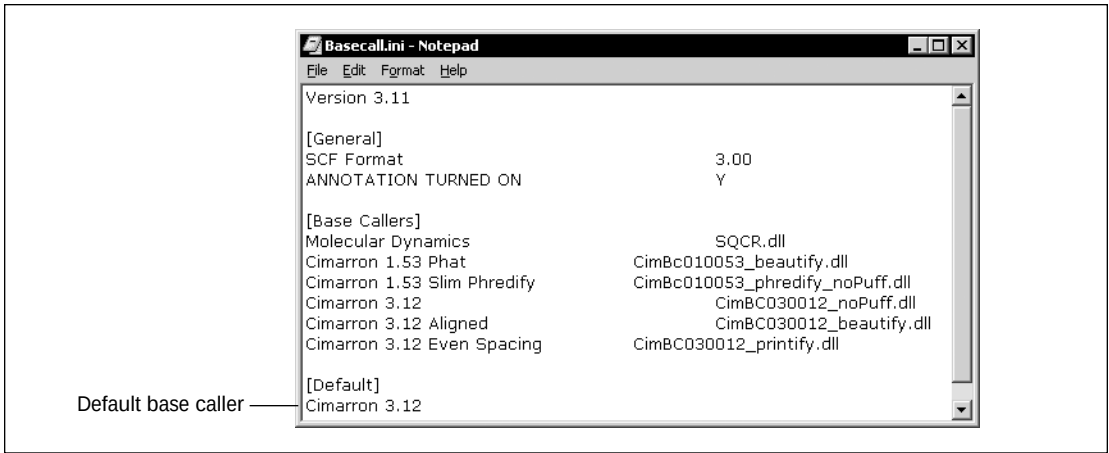


Figure 5-3. The Basecall.ini file.

Table 5-3. The sections and fields in the Basecall.ini file

Section and field name	Description	Setting type	Default
[General]			
SCF FORMAT	Determines the SCF file format version.	2.00 or 3.00	3.00
ANNOTATION TURNED ON	Determines whether the base caller produces annotation information, which is needed by the Sequencing Diagnostics software. Note: The base calling is slower if the setting is Y instead of N.	Y or N	Y
[Base Callers]			
	Lists the available base callers. The first column lists the base caller names, and the second column lists the corresponding file names of the base callers.	–	–
[Default]	Specifies the default base caller for the laboratory.	String*	Name of the latest base caller

* The text must match the name of one of the available base callers in the Base Callers section of the file.

5.7 Specifying which base caller to use

There are a variety of ways in which you can specify the base caller(s) to use. You can specify a—

- Default base caller for all plates
- Base caller(s) as part of the plate definition

5.7.1 Specifying a default base caller

If the plate definition for a given run does not specify a base caller, the software uses the default base caller. The software allows you to specify a default base caller in the Basecall.ini file (figure 5-3) and in the Automatic Base Calling window.

The Automatic Base Calling window

To open the Automatic Base Calling window (figure 5-2), on the **Configure** menu, click **Automatic Base Calling**. See the *MegaBACE 4000 Instrument User’s Guide* for details on the Automatic Base Calling window.

The default base caller specified in the Automatic Base Calling window overwrites the default base caller specified in the Basecall.ini file. When you restart the software, the default resets to the base caller specified in the Basecall.ini file.

Specifying the default base caller in the Basecall.ini file

The Base Callers section of the Basecall.ini file lists the names of the available base callers. The base callers are provided during installation. The Default section lists the name of the default base caller. To specify a default base caller, you type the base caller name in column one of the Default section (figure 5-3). The software recognizes only the base caller names listed in column one of the Base Callers section.

5.7.2 About specifying a base caller(s) in the plate definition

The Plate Setup window, a plate setup template (.tpl), and/or a plate setup data file (.psd) can be used in creating a plate definition. Each of these methods allows you to specify a base caller or base callers. If you want to specify—

- **Base caller(s) for multiple plates**—Use a plate setup template (.tpl) (chapter 6) and the Plate Setup window to create plate definitions.
- **Plate-specific and well-specific base callers**—Use a .psd file to create the plate definition. A .psd file allows you to define both a plate-level base caller(s) and well-level base callers for given wells on the plate (section 7.3).

If you use a combination of methods when creating plate definitions (for example, .tpl and .psd files) and the plate-level base caller is defined in more than one source, the software uses the plate-level base caller specified last.

Chapter 6 Creating templates for sequencing plate setup parameters

For each sequencing run, the Instrument Control Manager requires a plate definition and the instrument control parameters. This chapter describes plate definitions and how to create the plate setup templates for sequencing. The topics are—

- About creating the plate definitions (section 6.1)
- About plate setup templates for sequencing (section 6.2)
- Using the Plate Setup window to create a plate setup template (section 6.3)
- Creating a plate setup template in a text file (section 6.4)
- Using a template to suppress .rsd file creation for empty capillaries (section 6.5)
- Setting and clearing the default plate setup template (section 6.6)

For a description of the instrument control templates, see chapter 4.

6.1 About creating the plate definitions

A plate definition consists of the—

- **Plate ID**—The software uses the plate ID to name the run folder.
- **Plate setup parameters**—The electrophoresis parameters, chemistry parameters, sample names, comments and optional parameters for the plate.

6.1.1 Methods for creating the plate definitions

The Instrument Control Manager allows you to create the plate definitions in various ways. You can use—

- A plate setup data file (.psd) that an operator can import at the start of the Inject Samples and Run protocol. The .psd file can specify the plate ID and the plate setup parameters. In addition, the .psd file can include the well-specific attributes. For a description of the .psd files for sequencing, see chapter 7.
- The Plate Setup window (figure 6-1) to create the definitions before the run. The operator can select the precreated plate definitions at the start of the Inject Samples and Run protocol.

6.1.2 The methods for creating plate definitions in the Plate Setup window

You or an operator can create the plate definition before the run using the Plate Setup window (figure 6-1). To do this, you or the operator can—

- **Import a .psd file**—You can use a .psd file to enter the plate definition. When the .psd file name is entered in the Plate ID box of the Plate Setup window, the software imports the information from the .psd file and automatically creates the plate definition.
- **Select a plate setup template (.tpl)**—You can use a template to enter the plate setup parameters that are visible in the Plate Setup window for multiple plates (section 6.2). In the MegaBACE.ini file, you can specify a default template that appears automatically when an operator creates a new plate definition in the Plate Setup window. The software includes templates for various sequencing chemistries.
- **Manually enter or modify the setup parameters**—To manually enter or change parameters in the Plate Setup window, you must enable the edit mode in the MegaBACE.ini file. See section 2.3 for a description of the edit mode.

See the *MegaBACE 4000 Instrument User's Guide* for a description of how to create plate definitions in the Plate Setup window.

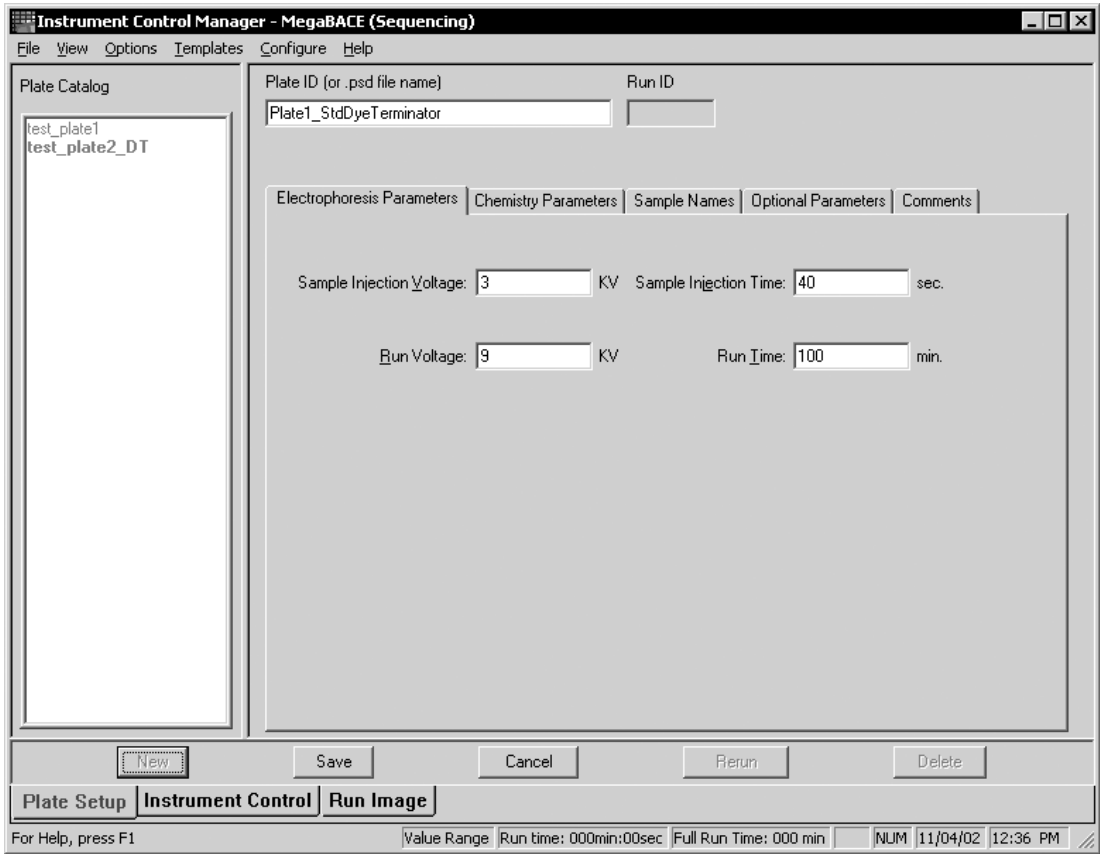


Figure 6-1. The Plate Setup window for sequencing.

6.2 About plate setup templates for sequencing

Important

The edit mode must be turned on to allow you to create, edit, select, or save plate setup templates (.tpl) in the Instrument Control Manager. See section 2.3 for information about enabling the edit mode.

You can simplify the process of creating plate definitions by using plate setup templates (.tpl files). In a plate setup template, you specify parameters that are common to multiple plates. You can create templates in the Plate Setup window (section 6.3), or you can create a .tpl file using any text editor or Excel.

See table 6-1 for a description of the plate setup parameters for sequencing.

Table 6-1. The recognized fields for sequencing plate setup templates

Field name	Value
The values for the following fields can be imported into the Plate Setup window or typed manually. When importing the values, the fields can be defined in either a .psd file or a plate setup template (.tpl).	
Note: If both a .psd file and a plate setup template are used to define the same parameter, the values provided in the file imported last overwrite the values for any identical fields.	
Electrophoresis Parameters tab	
Injection voltage	Integer or fraction (kV). The range is 1–15 kV.
Injection time	Integer (s). The range is 1–600 s.
Run voltage	Integer or fraction (kV). The range is 1–15 kV.
Run time	Integer (min). The range is 1–720 min.
Chemistry Parameters tab	
Chemistry	String. The name of a chemistry parameter set. Note: The name must match a name in the Chemistry.ini file (chapter 3) for the sequencing application. Otherwise, the software will not recognize the name.
Sample Names tab	
Well_ID (for example, A01) or Well_ID.NAME (.psd file only)	String. The sample name that you want to assign to each well location (ID) and use as the .rsd file name (default), for example, Sample1. Note: If you want to use names other than sample names as the file names, you must also use the Well_ID.FILE field in the .psd file.
Optional Parameters tab	
Note: The Run temperature and PMT voltages are typically defined as part of the instrument control parameters (chapter 4). The values can be defined here for a special run. If defined here, the values temporarily overwrite the same values defined in an .icp file or the Instrument Control window.	
Run temperature	Integer (°C). The range is 27–44 °C.
PMT1 voltage	Integer or fraction (kV). The range is 500–900 V.
PMT2 voltage	Integer or fraction (kV). The range is 500–900 V.
PMT3 voltage	Integer or fraction (kV). The range is 500–900 V.
PMT4 voltage	Integer or fraction (kV). The range is 500–900 V.

Table 6-1. The recognized fields for sequencing plate setup templates (continued)

Field name	Value
Base caller	String. The optional name of a plate-specific base caller that overwrites the default base caller specified in the Basecall.ini file or the Automatic Base Calling window. Note: You can designate multiple base callers in the .psd file or .tpl file for specified plates. You must separate each base caller name with a comma in the same value column for the BASE CALLER field.
Comments tab	
Comments	String, including spaces, representing the comments about the plate that you want to include.

6.3 Using the Plate Setup window to create a plate setup template

You can use the Plate Setup window (figure 6-1) to create a plate setup template. To do this—

1. In the Plate Setup window, click **New**.
2. In the following tabs, enter the parameters that you want to use for the template. See table 6-1 for a description of the parameters.
 - **Electrophoresis Parameters tab**—Type the values for the parameters in each of the boxes.
 - **Chemistry Parameters tab**—Select the name of the chemistry parameters from the **Chemistry Name** list.
 - **Sample Names tab**—By default, the software uses the well IDs to name the files (for example, A01.rsd). If you want to use sample names other than the well IDs, you can type the names in the Sample Names tab. You can also specify any empty wells for which you want to suppress the raw sample data file (.rsd) creation, see section 6.5.
 - **Optional Parameters tab**—(Optional) If you want to overwrite the Run Temperature and PMT Voltage values, which are provided in the instrument control parameters, you can type the values in this tab. You can also specify a base caller for the plate in this tab.
 - **Comments tab**—You can type any comments about the plate.

See the *MegaBACE 4000 Instrument User's Guide* for details on using the Plate Setup window.

- 3. (Optional) If you want to suppress the raw sample data file (.rsd) creation for empty capillaries, see section 6.5.
- 4. On the **Templates** menu, point to **Plate Setup Templates**, and click **Save Template**. The Save Template window appears.
- 5. Type the template name, and click **Save**. The template is saved in the ...\\MegaBACE\\Instrument Control\\Templates folder. The software adds the .tpl file extension.

6.4 Creating a plate setup template in a text file

Important When you create a plate setup template, you must type the .tpl file extension and save it in the ...\\MegaBACE\\Instrument Control\\Templates folder. Otherwise, the Instrument Control Manager will not find the file.

To create a plate setup template—

- 1. Open a new document or an existing .tpl file using a text editor (such as Notepad) or Excel.
- 2. Type the field names in the first column (field column), and type the values for each field in the second column (value column). Make sure the field column and value column are separated by a tab. See table 6-1 for a list of the recognized fields for the plate setup parameters. Figure 6-2 shows an example of a plate setup template file (.tpl).
- 3. Save the file in an appropriate format. You must type the .tpl file extension.

Note: If you use Excel, save the file as a tab-delimited file.

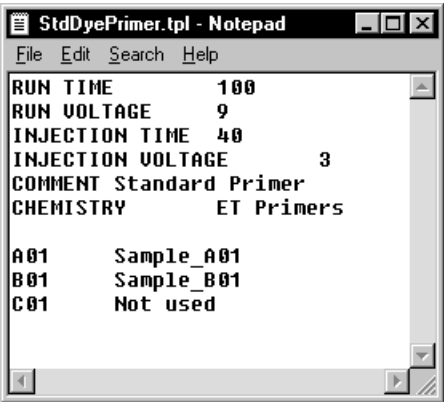


Figure 6-2. An example of a plate setup template file (.tpl) for sequencing.

6.5 Using a template to suppress .rsd file creation for empty capillaries

To specify an empty capillary, you designate the corresponding well as unused. The software does not create a raw sample data file (.rsd) for any wells that are specified as not used. In the Run Image window (figure 2-2), an unused well is indicated by a yellow arrow that appears above the capillary button.

Note: To specify bad capillaries, see section 2.6.

You can specify unused wells in either a—

- **Plate setup data file (.psd)**—You can use a .psd file to specify the unused wells if the unused wells occur for a single plate (section 7.6).
- **Plate setup template**—You can use a plate setup template if the unused wells occur for multiple plates.

You can specify unused wells in a plate setup template by typing **NOT USED** (case insensitive) as the sample name specified in either the Sample Names tab of the Plate Setup window or the .tpl file (section 6.4). To specify an unused well in the Sample Names tab of the Plate Setup window—

1. In the **Sample Names** tab of the Plate Setup window (figure 6-3), adjust the width of the column that includes the unused well. To increase the column width, click the **RIGHT ARROW** key in the Sample Names tab. To decrease the column width, click the **LEFT ARROW** key in the Sample Names tab.
2. Type **NOT USED** (case insensitive) as the sample name for the unused well.

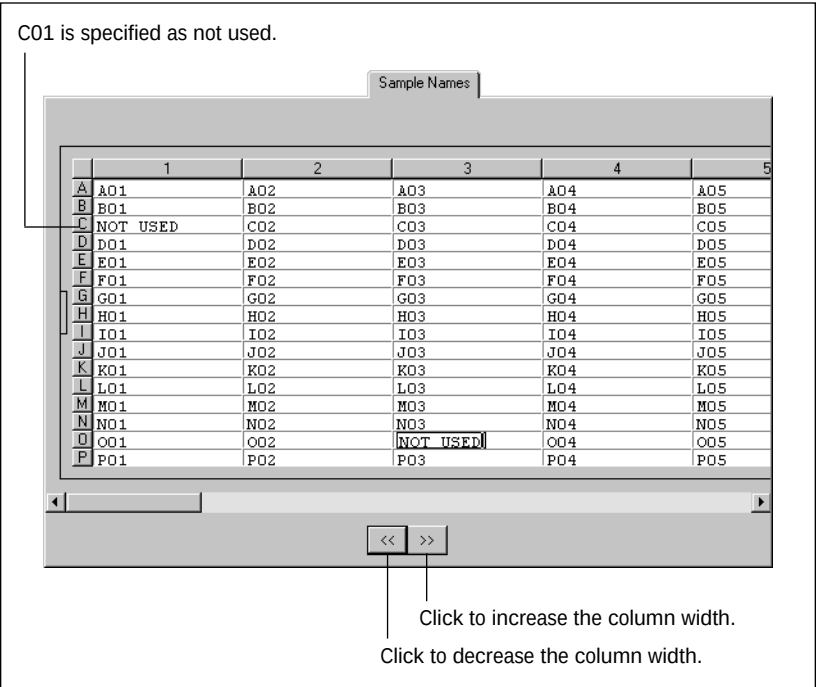


Figure 6-3. The Sample Names tab of the Plate Setup window.

6.6 Setting and clearing the default plate setup template

You can set a default plate setup template in the—

- **Sequencing section of the MegaBACE.ini file**—To set a default plate setup template for sequencing in the MegaBACE.ini file, see section 5.2.1.
- **Templates menu in the Plate Setup window**—The default template is active only until you close the Instrument Control Manager. When the software is restarted, the default plate setup template specified in the Sequencing section of the MegaBACE.ini file is active until you change it again.

6.6.1 Using the Templates menu to set a default template

To set the default plate setup template—

1. In the Instrument Control Manager, click the **Plate Setup** tab to display the Plate Setup window (figure 6-1).
2. On the **Templates** menu, point to **Plate Setup Template**, and click **Set Default**. The Set Default window appears.
3. Select the name of the template you want to use as the default, and click **Open**. The Instrument Control Manager designates the template you selected as the default template.

The next time you click New in the Plate Setup window, the values from the default template appear in the parameter boxes.

6.6.2 Clearing the default plate setup template

To clear the default plate setup template—

1. In the Instrument Control Manager, click the **Plate Setup** tab to display the Plate Setup window (figure 6-1).
2. On the **Templates** menu, point to **Plate Setup Templates**, and click **Clear Default**. The template you cleared is no longer the default template.

When you click New in the Plate Setup window, no values appear in the parameter boxes. You can select another default template (section 6.6.1), or use the Plate Setup window to create another template, save it, and then select it as the default template.

Important

When you restart the software, the default templates specified in the MegaBACE.ini file are active.

Chapter 7 Creating plate setup data files (.psd) for sequencing

This chapter describes the plate setup data files (.psd) for sequencing. The topics are—

- Using plate setup data files (.psd) for sequencing (section 7.1)
- Creating a plate setup data file (.psd) (section 7.2)
- Specifying a base caller in a .psd file (section 7.3)
- Annotating well-specific information in a .psd file (section 7.4)
- Designating the sample names and file names in the .psd file (section 7.5)
- Using a .psd file to suppress .rsd file creation for empty capillaries (section 7.6)

7.1 Using plate setup data files (.psd) for sequencing

You can define a plate ID, the plate setup parameters, and other plate- and well-level attributes in a .psd file so that an operator can automatically import the plate definition. With automatic import, the plate setup parameters are specified with less operator intervention.

The operator can import the .psd file either before the run in the Plate Setup window or during the run in the Select a Plate window (*MegaBACE 4000 Instrument User's Guide*). The parameters that you specify in the .psd file become part of the extended header of each raw sample data (.rsd) file created during the run.

7.1.1 About the plate ID

When an operator enters the text or bar code in the plate ID area of the Plate Setup window or the Select a Plate window (described in the *MegaBACE 4000 Instrument User's Guide*), the Instrument Control Manager looks for a .psd file with a file name that matches the text entered and automatically imports the plate setup parameters that are contained in the .psd file.

Part two Setting up the system for sequencing

An operator can enter a plate ID in the following ways:

- **Type the text** in the Plate ID box in the Plate Setup window or the Select a Plate window.
- **Scan a bar code** into the Plate ID box in the Plate Setup window or the Select a Plate window.

Important Regardless of the method used to enter the plate ID, the text entry or bar code (less a possible plate ID terminator) must match the name of the .psd file in order for a .psd file to be imported.

Figure 7-1 shows an example of a .psd file and the Plate Setup window. The .psd file name matches the plate ID typed in the Plate Setup window.

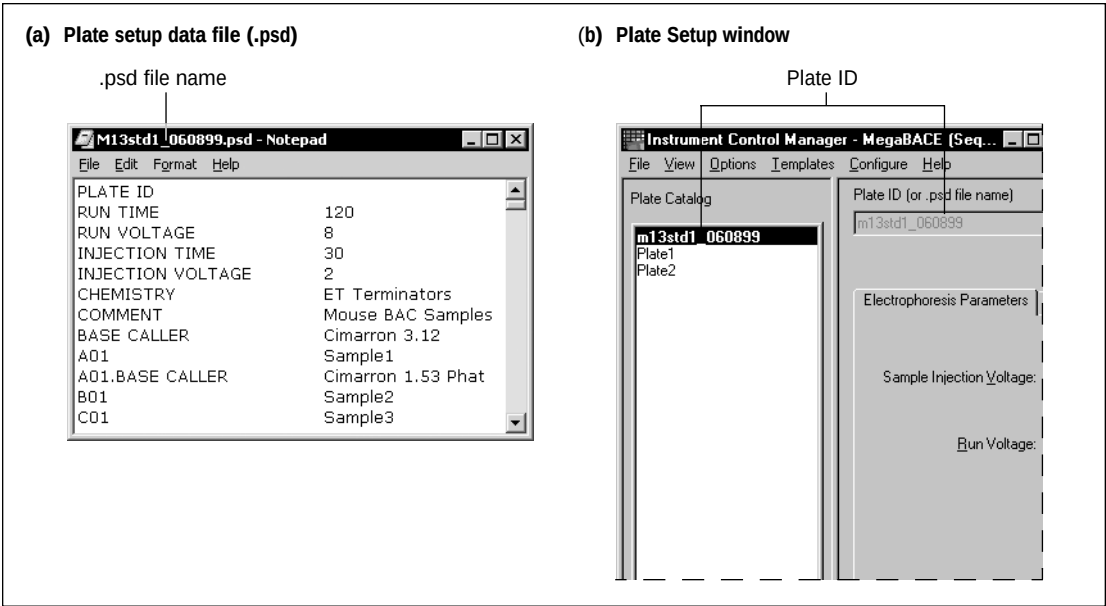


Figure 7-1. Using a .psd file: (a) Example of a .psd file. (b) The Plate Setup window with a plate ID that matches the .psd file name.

7.1.2 Using a plate ID that is different from the scanned ID

If the bar code is not what you want to use as the plate ID, you can use the .psd file to replace a bar code scanned into the Plate ID box in the Plate Setup window. To do this, you type a plate ID in the .psd file opposite the PLATE ID field (figure 7-2a).

The plate ID you type opposite the PLATE ID field in the .psd file replaces the bar code and becomes the real plate ID in the Plate Setup window (figure 7-2b).

For example, in figure 7-2, the—

- Plate ID scanned into the Plate Setup window is **10303385858**
- Plate setup data file name is **10303385858.psd**
- Plate ID opposite the PLATE ID field name inside the .psd file is **M13std40_092499** (figure 7-2a)
- Actual plate ID that appears in the Plate Setup window is **M13std40_092499** (figure 7-2b)

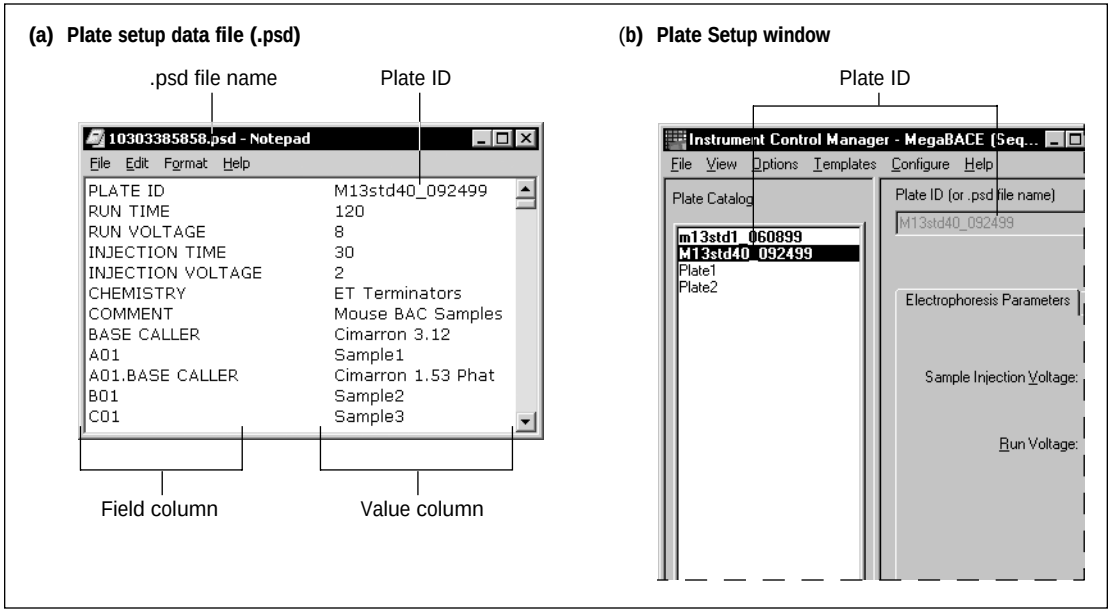


Figure 7-2. Replacing the scanned plate ID: (a) Example of a .psd file with a plate ID listed in the file. (b) The Plate Setup window with a plate ID that replaced the original bar code that was scanned.

7.2 Creating a plate setup data file (.psd)

Important

When you create a .psd file, you must type the .psd file extension and save it in the ...\\MegaBACE\\Instrument Control\\Psd folder (or the folder you specified as the Psd directory in the MegaBACE.ini file). Otherwise, the Instrument Control Manager will not find the file.

Regardless of the method the operator uses to enter the plate ID in the Plate Setup window or Select a Plate window, the text entry or bar code (less a possible plate ID terminator) must match the name of the .psd file in order for a .psd file to be imported.

To create a .psd file—

1. Open a new document or an existing .psd file using a text editor (such as Notepad) or Excel.
2. Type the field names in the first column (field column) and the values for each field name in the second column (value column). The order of the fields is arbitrary, and the field names are case-insensitive. Figure 7-3 shows an example of a .psd file created in Excel. Section 7.2.1 describes the fields that you can use in a .psd file.

Note: If you use spaces in the sample or file names, the Instrument Control Manager converts the spaces to underscore characters in the file name.

3. Save the file in an appropriate format. Make sure you use the .psd file extension.

Note: If you use Excel, save the file as a tab-delimited file (figure 7-3).

	A	B	C	D
1	PLATE ID	M13std40	092499	
2	RUN TIME	120		
3	RUN VOLTAGE	8		
4	INJECTION TIME	30		
5	INJECTION VOLTAGE	2		
6	CHEMISTRY	ET Terminators		
7	COMMENT	Mouse BAC Samples		
8	BASE CALLER	Cimarron 3.12		
9	A01	Sample1		
10	A01.BASE CALLER	Cimarron 1.53 Phat		
11	B01	NOT USED		
12	C01	Sample3		
13	D01	FB04_TSA_D01.x03		
14	E01	FB04_TSA_E01.x03		
15	F01	FB04_TSA_F01.x03		
16	G01	FB04_TSA_G01.x03		

Field column Value column

Figure 7-3. An example of a .psd file created in Excel.

7.2.1 The .psd field names for sequencing

The .psd file can include the following fields:

- All the fields that appear in the tabs of the Plate Setup window or in a plate setup template (.tpl). Table 6-1 lists the available fields.
Note: If both a .psd file and a plate setup template (.tpl) are used to define the same plate setup parameters, the values provided in the file imported last overwrite the values for any identical fields.
- Additional plate-specific fields and well-specific fields listed in table 7-1.

Table 7-1. Additional fields available for the sequencing .psd files

Field name	Value
Plate ID	String. If the plate ID entered here is different from the bar code of the .psd file name, the plate ID entered here becomes the real plate ID that appears in the Plate Setup window.
Well_ID.FILE (for example, A01.FILE)	String. The .rsd file name that you want to assign to a given well location. For example, A01 if you want to use the well location as the file name. This parameter overwrites the sample name as the file name.
Well_ID.BASE CALLER (for example, A01.BASE CALLER)	String. The name(s) of the base caller(s) that you want to use for a well that differ from the plate-level base caller(s). To specify multiple base callers, you separate each name with a comma.
User-defined fields	Fields that can be imported, but are not used by the software. These data are stored in the extended header of the files.

Tables 6-1 and 7-1 list the field names that the Instrument Control Manager recognizes, except for the user-defined fields. For additional plate information you want included in the output file (that will not appear in the Plate Setup window), you can use field names the Instrument Control Manager does not recognize (for example, Operator). You must also provide a name in the value column (for example, Mary Smith).

7.2.2 The .psd file value constraints

If you type a value that is invalid (for example, an alpha character when the software requires a number, or a value that is out of range), the Instrument Control Manager displays an error message when you import the .psd file. If this occurs, you must correct the value in the .psd file before it can be imported successfully. Alternatively, if the edit mode is on, you can enter the correct value in the Plate Setup window.

7.3 Specifying a base caller in a .psd file

You can specify a base caller or multiple base callers for a plate and for wells on a plate in the .psd file. The base caller(s) specified in the .psd file overwrite the base caller(s) specified in the default plate setup template (.tpl) (section 6.4) and the default base caller specified in the Basecall.ini file (section 5.6).

7.3.1 Designating a base caller for the plate

To designate a plate-specific base caller in a .psd file, you use the following field:

Field name	Value
BASE CALLER	The base caller name(s)

To designate a plate-level base caller, you can use the BASE CALLER field in either the .psd file or the .tpl file. To designate multiple plate-level base callers, separate each base caller name with a comma in the value column of the BASE CALLER field.

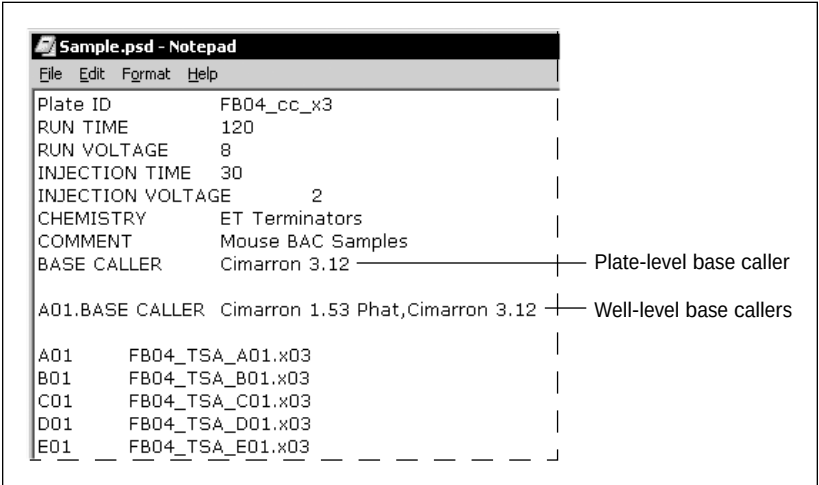


Figure 7-4. Example of a .psd file with entries to designate plate- and well-specific base callers.

7.3.2 Specifying a base caller for a selected well on a plate

To specify one or more well-level base callers that differ from the plate-level base caller, you use the following field in the .psd file. If you specify multiple base callers, separate each base caller name with a comma (figure 7-4).

Field name	Value
Well_ID.BASE CALLER	The base caller name(s)

7.4 Annotating well-specific information in a .psd file

The values for the field names that begin with a well ID (for example, A01) are well specific. For example, you can specify the following types of information for a given well:

- Sample name, which becomes the .rsd file name by default
- Base callers (section 7.3)
- Explicit file name, if you need the .rsd file name to be different from the sample name (section 7.5)
- User-defined information

If you want additional well-specific information for your records, the Instrument Control Manager allows you to include a dot (.) character after the well ID followed by an additional name to annotate the well ID. The following and figure 7-2 show examples of what you might enter in the .psd file for annotation of well-specific information:

Field name	Value
A01.BASE CALLER	Cimarron 1.53 Phat,Cimarron 3.12
A01.SAMPLE TYPE	Tissue

Note: To designate multiple well-specific base callers, separate each base caller name with a comma in the .psd file.

The Instrument Control Manager writes the additional well-specific information only into the .rsd file header for the corresponding well. The well-specific annotation does not appear in the Plate Setup window.

7.5 Designating the sample names and file names in the .psd file

If no file names are specified, but the sample names are provided for a plate, the software uses the sample names as the file names, for example, sample1.rsd. If no sample or file names are provided for a plate, the Instrument Control Manager uses the well location (well IDs) as the sample file names, for example, A01.rsd.

You can provide sample names and designate explicit file names that are different from the sample names. To do this, you use two fields—Well_ID.FILE and Well_ID or Well_ID.NAME. The following and figure 7-5 show examples of what you might enter in the .psd file to designate a file name that is different from the sample name.

Field name	Value	Comment
A01.NAME	14660	Sample name
A01.FILE	A01	File name

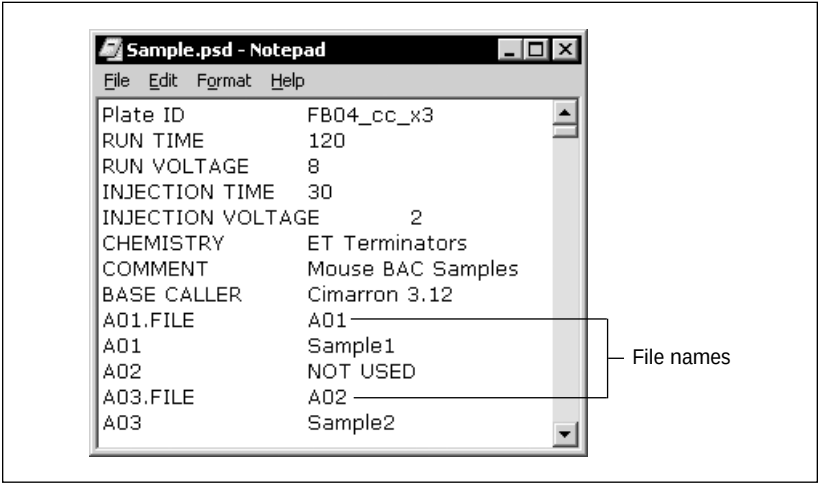


Figure 7-5. Example of a .psd file with entries to designate the well ID as the file names.

7.6 Using a .psd file to suppress .rsd file creation for empty capillaries

You can indicate unused wells by typing **NOT USED** (case-insensitive) as the sample name specified in the .psd file (figure 7-5). For every unused well, a yellow arrow appears above the capillary button in the Run Image window (figure 2-2). The software does not create an .rsd file for an unused well. For details on suppressing .rsd files for defective capillaries, see section 2.6.

Part three

Setting up the system for STR genotyping

Chapter 8 Specifying preferences for STR genotyping

This chapter describes the fields that have been added to the system initialization files to support STR genotyping. The topics are—

- The STR genotyping fields in the MegaBACE.ini file (section 8.1)
- Setting STR genotyping as the default application (section 8.2)
- Setting default parameter templates for STR genotyping (section 8.3)

Make sure you configure the General section of the MegaBACE.ini file (chapter 2) for your laboratory. To create additional chemistry parameters, see chapter 3.

8.1 The STR genotyping fields in the MegaBACE.ini file

The MegaBACE.ini file determines the settings available in the Instrument Control Manager when you open the software. The file is divided into sections: General, Sequencing, Genotyping (STR genotyping), and SNP (SNP genotyping). Figure 8-1 shows the MegaBACE.ini file with the Genotyping fields highlighted.

Table 8-1 describes the fields in the Genotyping section of the MegaBACE.ini file. See chapter 2 for a description of the other fields in the MegaBACE.ini file.

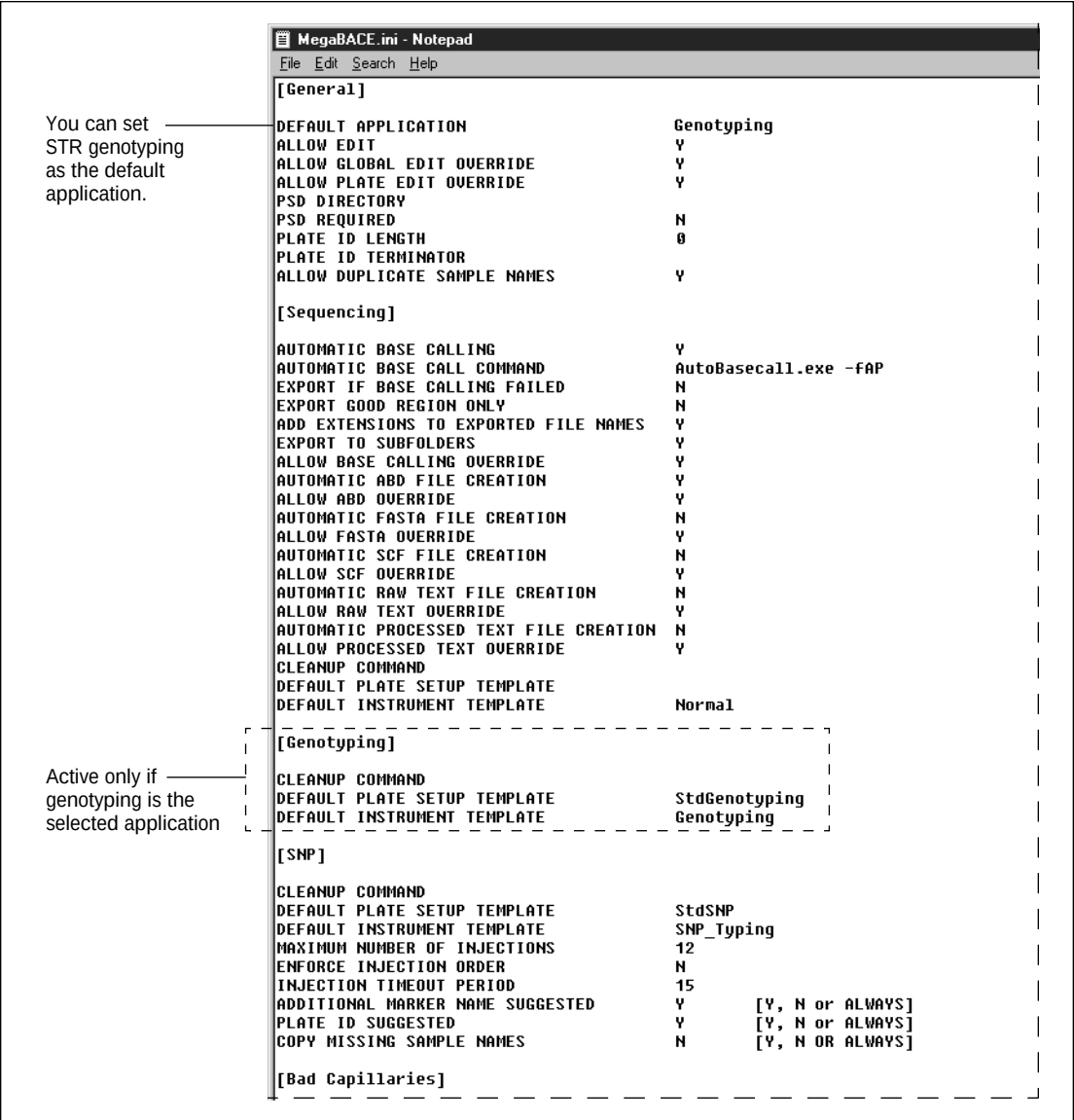


Figure 8-1. The STR genotyping fields in the MegaBACE.ini file. Note that the file contains a section for each application, which allows you to specify a default plate setup and instrument template.

Table 8-1. The STR genotyping fields in the MegaBACE.ini file

Field name	Description	Setting type	Default
Cleanup command	Allows the administrator to specify a program or batch file to be executed at the end of the run. This can be useful for transferring the raw data files to another workstation for analysis or archiving. See appendix C for a more detailed description.	String consisting of the path to the program or batch file (full or relative path to the instrument control folder)	none
Default plate setup template	Specifies a default plate setup template.	String	StdGenotyping
Default instrument template	Specifies a default instrument control template.	String	Genotyping

8.2 Setting STR genotyping as the default application

The operator can change the selected application in the Instrument Control Manager by choosing the Applications command on the View menu. However, when you restart the Instrument Control Manager, the software automatically resets the selected application to the default specified in the MegaBACE.ini file.

To set STR genotyping as the default application, you use the following field and value in the MegaBACE.ini file:

Field in MegaBACE.ini file	Value
DEFAULT APPLICATION	Genotyping

8.3 Setting default parameter templates for STR genotyping

The Instrument Control Manager allows you to use templates to define the plate setup parameters and the instrument control parameters for multiple runs. You can specify default templates in the Genotyping section of the MegaBACE.ini file. The software includes the following templates for the STR genotyping application, which are specified as the default templates:

- **Plate setup template**—StdGenotyping.tpl
- **Instrument control template**—Genotyping.icp



If you enable the edit mode in the MegaBACE.ini file (section 2.3), the operator can modify the values for a given run or choose a different template in the Instrument Control Manager. However, when the Instrument Control Manager is restarted, the software automatically resets the selected templates to the defaults specified in the MegaBACE.ini file.

You can use the templates provided, or modify the existing templates to create new templates. To create a new plate setup template for STR genotyping, see chapter 9. To create a new instrument control template, see chapter 4. To set new STR genotyping defaults in the MegaBACE.ini file, see the following sections.

8.3.1 Setting the default plate setup template for STR genotyping

If you specify a default template in the MegaBACE.ini file, the values from the plate setup template appear automatically in the Plate Setup window. As figure 8-1 and table 8-1 show, the software has been configured to use the StdGenotyping.tpl as the default template for STR genotyping. If you need to specify a different default template for the STR genotyping application, you can use the following field in the Genotyping section of the MegaBACE.ini file:

Field in MegaBACE.ini file	Value
DEFAULT PLATE SETUP TEMPLATE	Type the file name of the .tpl file, without the .tpl file extension.

If the edit mode is enabled, the operator can select a different default plate setup template from the Templates menu in the Plate Setup window.

Note: If the operator uses a plate setup data file (.psd) to define the plate for a run, the values in the .psd file overwrite any duplicate parameters specified in a .tpl file.

8.3.2 Setting the default instrument control template for STR genotyping

If you specify a default template in the MegaBACE.ini file, the values from the instrument control template appear automatically in the Instrument Control window when the operator performs a run. As figure 8-1 and table 8-1 show, the software has been configured to use the Genotyping.icp template as the default for STR genotyping. If you need to specify a different default template for the STR genotyping application, you can use the following field in the Genotyping section of the MegaBACE.ini file:



Field in MegaBACE.ini file	Value
DEFAULT INSTRUMENT TEMPLATE	Type the file name of the .icp file, without the .icp file extension.

If the edit mode is enabled, the operator can select a different instrument template from the Templates menu in the Instrument Control window.

Chapter 9 Creating templates for STR genotyping plate setup parameters

For each STR genotyping run, the Instrument Control Manager requires a plate definition and the instrument control parameters. This chapter describes plate definitions and how to create the plate setup templates for STR genotyping. The topics are—

- About creating the plate definitions (section 9.1)
- About plate setup templates for STR genotyping (section 9.2)
- Using the Plate Setup window to create a template (section 9.3)
- Creating a plate setup template in a text file (section 9.4)
- Using a template to suppress .rsd file creation for empty capillaries (section 9.5)
- Setting and clearing the default plate setup template (section 9.6)

For a description of the instrument control templates, see chapter 4.

9.1 About creating the plate definitions

A plate definition consists of the—

- **Plate set ID**—The software uses the plate ID to name the run folder.
- **Plate setup parameters**—The electrophoresis parameters, chemistry parameters, sample names, comments and optional parameters for the plate.

9.1.1 Methods for creating the plate definitions

The Instrument Control Manager allows you to create the plate definitions in various ways. You can use—

- A plate setup data file (.psd) that an operator can import at the start of the Inject Samples and Run protocol. The .psd file can specify the plate ID and the plate setup parameters. In addition, the .psd file can include the well-specific attributes. For a description of the .psd files for STR genotyping, see chapter 10.
- The Plate Setup window (figure 9-1) to create the definitions before the run. The operator can select the precreated plate definitions at the start of the Inject Samples and Run protocol.

9.1.2 The methods for creating plate definitions in the Plate Setup window

You or an operator can create the plate definition before the run using the Plate Setup window (figure 9-1). To do this, you or the operator can—

- **Import a .psd file**—You can use a .psd file to enter the plate definition specific to a single plate. The .psd file can specify plate attributes in addition to the parameters that are visible in the Plate Setup window. When the .psd file name is entered in the Plate ID box of the Plate Setup window, the software imports the information from the .psd file and automatically creates the plate definition.
- **Select a plate setup template (.tpl)**—You can use a template to enter the plate setup parameters for multiple plates (section 9.2). The template can specify only the parameters that are visible in the Plate Setup window. In the MegaBACE.ini file, you can specify a default template that appears automatically when an operator creates a new plate definition in the Plate Setup window. The software includes a default template for the STR genotyping application (StdGenotyping.tpl).
- **Manually enter or modify the setup parameters**—To manually enter or change parameters in the Plate Setup window, you must enable the edit mode in the MegaBACE.ini file. See section 2.3 for a description of the edit mode.

See the *MegaBACE 4000 Instrument User's Guide* for a description of how to create plate definitions in the Plate Setup window.

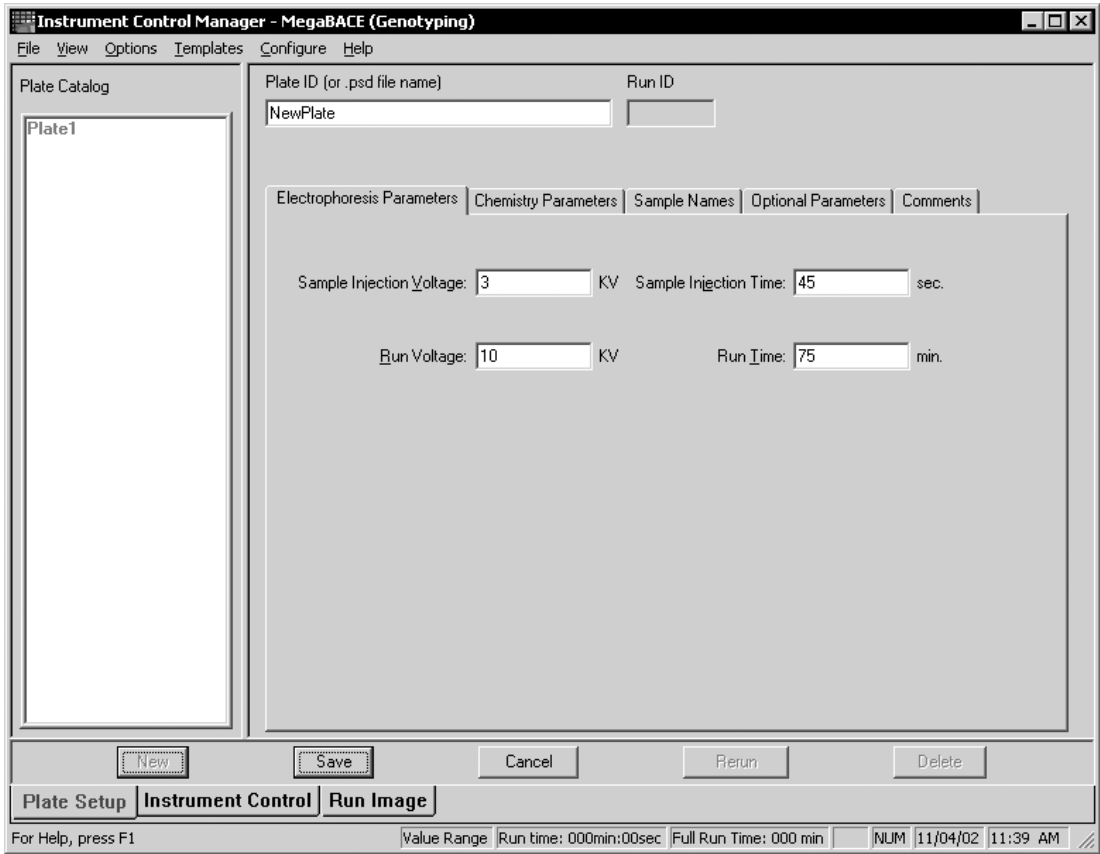


Figure 9-1. The Plate Setup window for the STR genotyping application.

9.2 About plate setup templates for STR genotyping

Important

The edit mode must be turned on to allow you to create, edit, select, or save plate setup templates (.tpl) in the Instrument Control Manager. See section 2.3 for information about enabling the edit mode.

You can simplify the process of creating plate definitions by using plate setup templates (.tpl). The templates allow you to save the frequently used plate setup parameters so that an operator does not have to re-enter the required parameters for every plate. In a plate setup template, you specify parameters that are common to multiple plates. You can create templates in the Plate Setup window (section 9.3), or you can create a .tpl file using any text editor or Excel (section 9.4).

See table 9-1 for a list of the recognized fields for the plate setup parameters.

Table 9-1. The recognized fields for the STR genotyping plate setup templates

Field name	Value
The values for the following fields can be imported into the Plate Setup window or typed manually. When importing the values, the fields can be defined in either a .psd file or a plate setup template (.tpl).	
Note: If both a .psd file and template are used to define the same parameter, the values in the file imported last overwrite the values for any identical fields.	
Electrophoresis Parameters tab	
Injection voltage	Integer or fraction (kV). The range is 1–15 kV.
Injection time	Integer (s). The range is 1–600 s.
Run voltage	Integer or fraction (kV). The range is 1–15 kV.
Run time	Integer (min). The range is 1–720 min.
Chemistry Parameters tab	
Chemistry	String. The name of a chemistry parameter set. Note: The name must match a name in the Chemistry.ini file (chapter 3) for the STR genotyping application. Otherwise, the software will not recognize the name.
Sample Names tab	
Well_ID (for example, A01) or Well_ID.NAME (.psd file only)	String. The sample name that you want to assign to each well location (ID) and use as the .rsd file name (default), for example, Sample1. Note: If you want to use names other than sample names as the file names, you must also use the Well_ID.FILE field in the .psd file.
Optional Parameters tab	
Note: The Run temperature and PMT voltages are typically defined as part of the instrument control parameters (chapter 4). The values can be defined here for a special run. If defined here, the values temporarily overwrite the same values defined in an .icp file or the Instrument Control window.	
Run temperature	Integer (°C). The range is 27–44 °C.
PMT1 voltage	Integer or fraction (kV). The range is 500–900 V.
PMT2 voltage	Integer or fraction (kV). The range is 500–900 V.
PMT3 voltage	Integer or fraction (kV). The range is 500–900 V.
PMT4 voltage	Integer or fraction (kV). The range is 500–900 V.
Comments tab	
Comments	String, including spaces, representing the comments about the plate that you want to include.

9.3 Using the Plate Setup window to create a template

You can use the Plate Setup window (figure 9-1) to create a plate setup template. To do this—

1. In the Plate Setup window, click **New**.
2. In the following tabs, enter the parameters that you want to use for the template. See table 9-1 for a description of the parameters.
 - **Electrophoresis Parameters tab**—Type the values for the parameters in each of the boxes.
 - **Chemistry Parameters tab**—Select the name of the chemistry parameters from the **Chemistry Name** list.
 - **Sample Names tab**—By default, the software uses the well IDs to name the files (for example, A01.rsd). If you want to use sample names other than the well IDs, you can type the names in the Sample Names tab.

If you want to use the Sample Names tab to specify empty wells for which you want to suppress the raw sample data file (.rsd) creation, see section 9.5.

- **Optional Parameters tab**—(Optional) If you want to overwrite the Run Temperature and PMT Voltage values, which are provided in the instrument control parameters, you can type the values in this tab.
- **Comments tab**—You can type any comments about the plate.

See the *MegaBACE 4000 Instrument User's Guide* for details on using the Plate Setup window.

3. On the **Templates** menu, point to **Plate Setup Templates**, and click **Save Template**. The Save Template window appears.
4. Type the template name, and click **Save**. The template is saved in the ...\\MegaBACE\\Instrument Control\\Templates folder. The software adds the .tpl file extension.

Important

9.4 Creating a plate setup template in a text file

When you create a plate setup template, you must type the .tpl file extension and save it in the ...MegaBACE\Instrument Control\Templates folder. Otherwise, the Instrument Control Manager will not find the file.

To create a plate setup template—

1. Open a new document or an existing .tpl file using a text editor (such as Notepad) or Excel.
2. Type the field names in the first column (field column), and type the values for each field in the second column (value column). Make sure the field column and value column are separated by a tab. See table 9-1 for a list of the recognized fields for the plate setup parameters. Figure 9-2 shows an example of a plate setup template file (.tpl).
3. Save the file in an appropriate format. You must type the .tpl file extension.

Note: If you use Excel, save the file as a tab-delimited file.

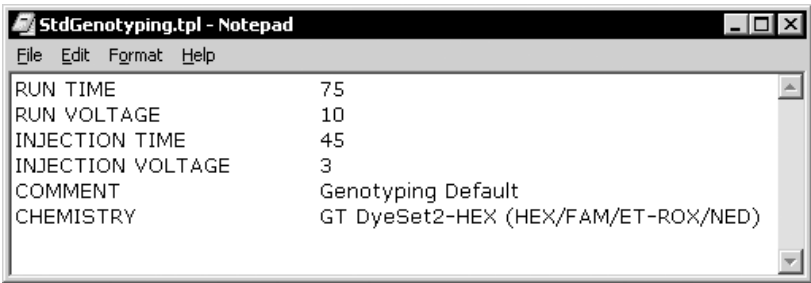


Figure 9-2. An example of a plate setup template file (.tpl) for STR genotyping.

9.5 Using a template to suppress .rsd file creation for empty capillaries

To specify an empty capillary, you designate the corresponding well as unused. The software does not create a raw sample data file (.rsd) for any wells that are specified as not used. In the Run Image window (figure 2-2), an unused well is indicated by a yellow arrow that appears above the capillary button.

Note: To specify bad capillaries, see section 2.6.

You can specify unused wells in either a—

- **Plate setup data file (.psd)**—You can use a .psd file to specify the unused wells if the unused wells occur for a single plate (section 7.6).

- **Plate setup template**—You can use a plate setup template if the unused wells occur for multiple plates.

You can specify unused wells in a plate setup template by typing **NOT USED** (case-insensitive) as the sample name specified in either the Sample Names tab of the Plate Setup window or the .tpl file (section 9.4).

To specify an unused well in the Sample Names tab of the Plate Setup window—

1. In the **Sample Names** tab of the Plate Setup window (figure 9-3), adjust the width of the column that includes the unused well. To increase the column width, click the **RIGHT ARROW** key in the Sample Names tab. To decrease the column width, click the **LEFT ARROW** key in the Sample Names tab.
2. Type **NOT USED** (case-insensitive) as the sample name for the unused well.

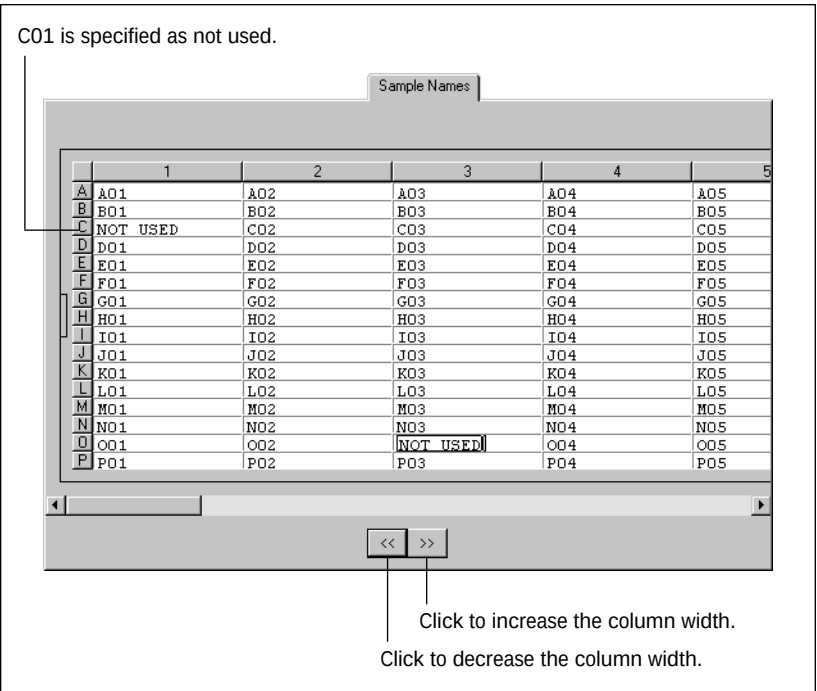


Figure 9-3. The Sample Names tab of the Plate Setup window.

9.6 Setting and clearing the default plate setup template

You can set a default plate setup template in the—

- **Genotyping section of the MegaBACE.ini file**—To set a default plate setup template for STR genotyping in the MegaBACE.ini file, see section 8.3.1.
- **Templates menu in the Plate Setup window**—The default template is active only until you close the Instrument Control Manager. When the software is restarted, the default plate setup template specified in the Genotyping section of the MegaBACE.ini file is active until you change it again.

9.6.1 Using the Templates menu to set a default template

To set the default plate setup template—

1. In the Instrument Control Manager, click the **Plate Setup** tab to display the Plate Setup window (figure 9-1).
2. On the **Templates** menu, point to **Plate Setup Template**, and click **Set Default**. The Set Default window appears.
3. Select the name of the template you want to use as the default, and click **Open**. The Instrument Control Manager designates the template you selected as the default template.

The next time you click New in the Plate Setup window, the values from the default template appear in the parameter boxes.

9.6.2 Clearing the default plate setup template

To clear the default plate setup template—

1. In the Instrument Control Manager, click the **Plate Setup** tab to display the Plate Setup window (figure 9-1).
2. On the **Templates** menu, point to **Plate Setup Templates**, and click **Clear Default**. The template you cleared is no longer the default template.

When you click New in the Plate Setup window, no values appear in the parameter boxes. You can select another default template (section 9.6.1), or use the Plate Setup window to create another template, save it, and then select it as the default template.

Important

When you restart the software, the default templates specified in the MegaBACE.ini file are active.

Chapter 10 Creating plate setup data files (.psd) for STR genotyping

This chapter describes the plate setup data files (.psd) for STR genotyping. The topics are—

- Using plate setup data files (.psd) for STR genotyping (section 10.1)
- Creating a plate setup data file (.psd) (section 10.2)
- Annotating well-specific information in a .psd file (section 10.3)
- Designating the sample names and file names in the .psd file (section 10.4)
- Using a .psd file to suppress .rsd file creation for empty capillaries (section 10.5)

10.1 Using plate setup data files (.psd) for STR genotyping

You can define a plate ID, the plate setup parameters, and other plate- and well-level attributes in a .psd file so that an operator can automatically import the plate definition. With automatic import, the plate setup parameters are specified with less operator intervention.

The operator can import the .psd file either before the run in the Plate Setup window or during the run in the Select a Plate window (*MegaBACE 4000 Instrument User's Guide*). The parameters that you specify in the .psd file become part of the extended header of each raw sample data (.rsd) file created during the run.

10.1.1 About the plate ID

When an operator enters the text or bar code in the plate ID area of the Plate Setup window or the Select a Plate window (described in the *MegaBACE 4000 Instrument User's Guide*), the Instrument Control Manager looks for a .psd file with a file name that matches the text entered and automatically imports the plate setup parameters that are contained in the .psd file.

Part three Setting up the system for STR genotyping

An operator can enter a plate ID in the following ways:

- **Type the text** in the Plate ID box in the Plate Setup window or the Select a Plate window.
- **Scan a bar code** into the Plate ID box in the Plate Setup window or the Select a Plate window.

Important Regardless of the method used to enter the plate ID, the text entry or bar code (less a possible plate ID terminator) must match the name of the .psd file in order for a .psd file to be imported.

Figure 10-1 shows an example of a .psd file and the Plate Setup window. The .psd file name matches the plate ID typed in the Plate Setup window.

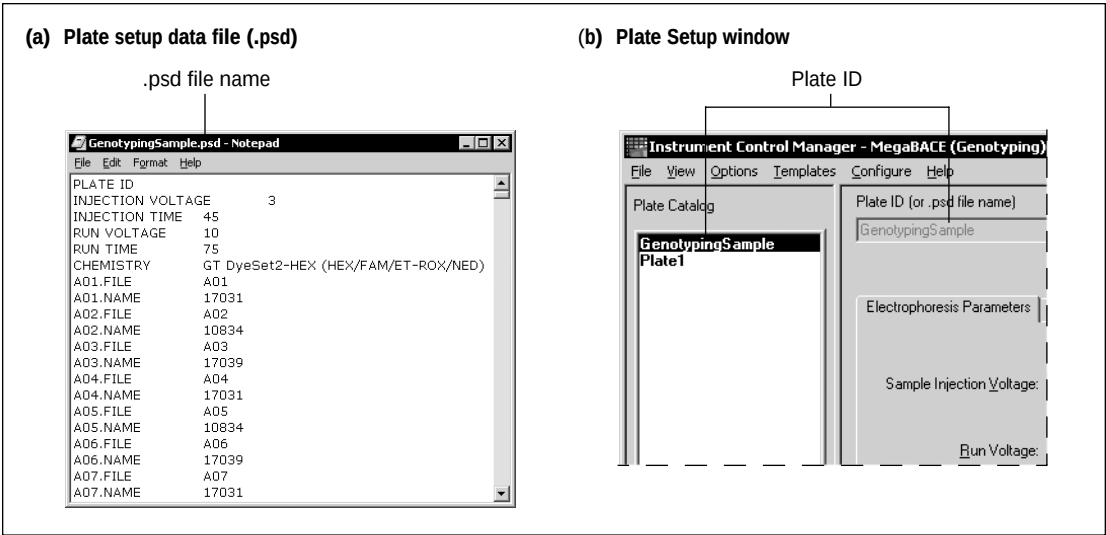


Figure 10-1. Using a .psd file: (a) Example of a .psd file. (b) The Plate Setup window with a plate ID that matches the .psd file name.

10.1.2 Using a plate ID that is different from the scanned ID

If the bar code is not what you want to use as the plate ID, you can use the .psd file to replace a bar code scanned into the Plate ID box in the Plate Setup window. To do this, you type a plate ID in the .psd file opposite the PLATE ID field (figure 10-2a).

The plate ID you type opposite the PLATE ID field in the .psd file replaces the bar code and becomes the real plate ID in the Plate Setup window (figure 10-2b).

For example, in figure 10-2, the—

- Plate ID scanned into the Plate Setup window is **GenotypingSample**
- Plate setup data file name is **GenotypingSample.psd**
- Plate ID opposite the PLATE ID field name inside the .psd file is **Genotyping_Demo_Panel6** (figure 10-2a)
- Actual plate ID that appears in the Plate Setup window is **Genotyping_Demo_Panel6** (figure 10-2b)

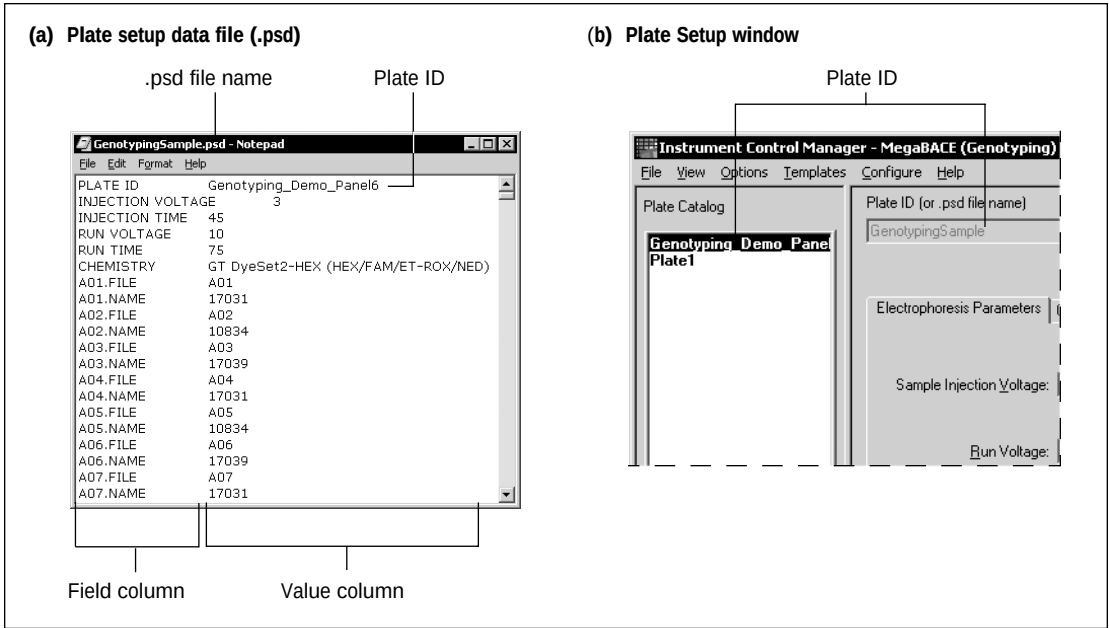


Figure 10-2. Replacing the scanned plate ID: (a) Example of a .psd file with a plate ID listed in the file.
(b) The Plate Setup window with a plate ID that replaced the original bar code that was scanned.

10.2 Creating a plate setup data file (.psd)

Important

When you create a .psd file, you must type the .psd file extension and save it in the ...\\MegaBACE\\Instrument Control\\Psd folder (or the folder you specified as the Psd directory in the MegaBACE.ini file). Otherwise, the Instrument Control Manager will not find the file.

Regardless of the method the operator uses to enter the plate ID in the Plate Setup window or Select a Plate window, the text entry or bar code (less a possible plate ID terminator) must match the name of the .psd file in order for a .psd file to be imported.


To create a .psd file—

1. Open a new document or an existing .psd file using a text editor (such as Notepad) or Excel.
2. Type the field names in the first column (field column) and the values for each field name in the second column (value column). The order of the fields is arbitrary, and the field names are case-insensitive. Figure 10-3 shows an example of a .psd file created in Excel. Section 10.2.1 describes the fields that you can use in a .psd file.

Note: If you use spaces in the sample or file names, the Instrument Control Manager converts the spaces to underscore characters in the file name.

3. Save the file in an appropriate format. You must type the .psd file extension.

Note: If you use Excel, save the file as a tab-delimited file (figure 10-3).

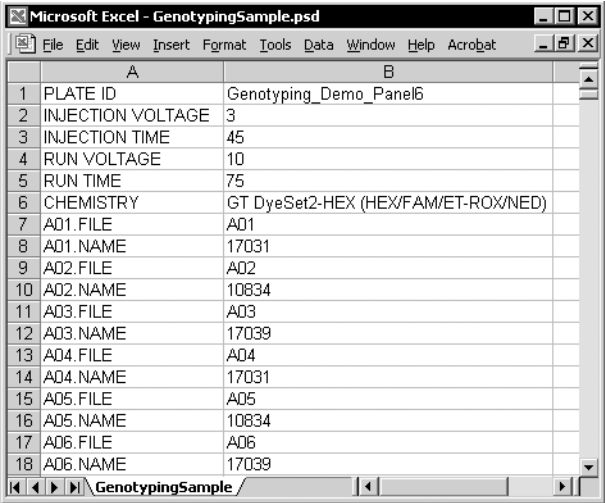

	A	B
1	PLATE ID	Genotyping_Demo_Panel6
2	INJECTION VOLTAGE	3
3	INJECTION TIME	45
4	RUN VOLTAGE	10
5	RUN TIME	75
6	CHEMISTRY	GT DyeSet2-HEX (HEX/FAM/ET-ROX/NED)
7	A01.FILE	A01
8	A01.NAME	17031
9	A02.FILE	A02
10	A02.NAME	10834
11	A03.FILE	A03
12	A03.NAME	17039
13	A04.FILE	A04
14	A04.NAME	17031
15	A05.FILE	A05
16	A05.NAME	10834
17	A06.FILE	A06
18	A06.NAME	17039

Field column Value column

Figure 10-3. An example of a .psd file created in Excel.

10.2.1 The .psd field names

The .psd file can include the following fields:

- All the fields that appear in the tabs or Plate Setup window or in a plate setup template (.tpl). Table 9-1 lists the available fields.

Note: If both a .psd file and a plate setup template (.tpl) are used to define the same plate setup parameters, the values provided in the file imported last overwrite the values for any identical fields.

- The plate- and well-specific fields listed in table 10-1.

Table 10-1. Additional fields available for the STR genotyping .psd files

Field name	Value
Plate ID	The string for the plate ID. If the plate ID entered here is different from the bar code of the .psd file name, the plate ID entered here becomes the real plate ID that appears in the Plate Setup window.
Well_ID.FILE (for example, A01.FILE)	The .rsd file name that you want to assign for a given well location. For example, A01 if you want to use the well location as the file name. This parameter overwrites the sample name as the file name.
User-defined fields	Fields that can be imported, but are not used by the software. These data are stored in the extended header of the files.

Tables 9-1 and 10-1 list the field names that the Instrument Control Manager recognizes, except for the user-defined fields. For additional plate information you want included in the output file (that will not appear in the Plate Setup window), you can use field names the Instrument Control Manager does not recognize (for example, Operator). You must also provide a name in the value column (for example, Mary Smith).

10.2.2 The .psd file value constraints

If you type a value that is invalid (for example, an alpha character when the software requires a number, or a value that is out of range), the Instrument Control Manager displays an error message when you import the .psd file. If this occurs, you must correct the value in the .psd file before it can be imported successfully. Alternatively, if the edit mode is on, you can enter the correct value in the Plate Setup window.

10.3 Annotating well-specific information in a .psd file

The values for the field names that begin with a well ID (for example, A01) are well specific. If you want additional well-specific information for your records, the Instrument Control Manager allows you to include a dot (.) character after the well ID followed by an additional name to annotate the well ID. Figure 10-3 shows examples of what you might enter in the .psd file for annotation of well-specific information.

The Instrument Control Manager writes the additional well-specific information only into the .rsd file header for the corresponding well. The well-specific annotation does not appear in the Plate Setup window.

10.4 Designating the sample names and file names in the .psd file

If no file names are specified, but the sample names are provided for a plate, the software uses the sample names as the file names, for example, sample1.rsd. If no sample or file names are provided for a plate, the Instrument Control Manager uses the well location (well IDs) as the sample file names, for example, A01.rsd.

You can provide sample names and designate explicit file names that are different from the sample names. To do this, you use two fields—Well_ID.FILE and Well_ID.NAME. The following and figure 10-4 show examples of what you might enter in the .psd file to designate a file name that is different from the sample name.

Field name	Value	Comment
A01.NAME	14660	Sample name
A01.FILE	A01	File name

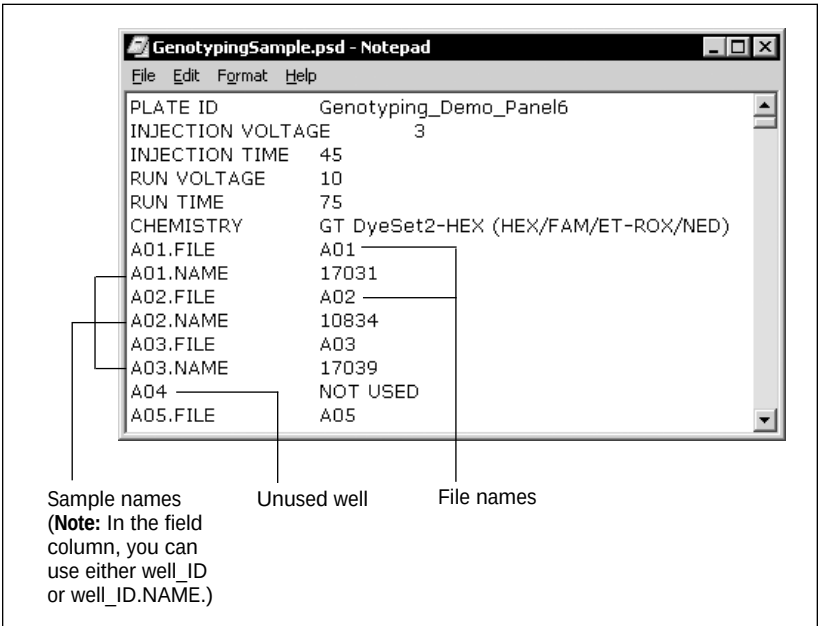


Figure 10-4. Example of a .psd file with entries to designate the well ID as the file names.

10.5 Using a .psd file to suppress .rsd file creation for empty capillaries

You can indicate unused wells by typing **NOT USED** (case-insensitive) as the sample name specified in the .psd file (figure 10-4). For every unused well, a yellow arrow appears above the capillary button in the Run Image window (figure 2-2). The software does not create an .rsd file for an unused well. For details on suppressing .rsd files for defective capillaries, see section 2.6.

Part four

Setting up the system for SNP genotyping

Chapter 11 Specifying preferences for SNP genotyping

This chapter describes the fields in the MegaBACE.ini file that you can use to configure the software for SNP genotyping. The topics are—

- The SNP genotyping fields in the MegaBACE.ini file (section 11.1)
- Setting SNP genotyping as the default application (section 11.2)
- Setting default parameter templates for SNP genotyping (section 11.3)
- Specifying limits for the injections (section 11.4)
- Enforcing the order in which an operator injects the plates (section 11.5)
- Prompting the operator for plate-specific information (section 11.6)

Make sure you configure the General section of the MegaBACE.ini file for your laboratory (chapter 2). To create additional chemistry parameters, see chapter 3.

11.1 The SNP genotyping fields in the MegaBACE.ini file

The MegaBACE.ini file determines the settings available in the Instrument Control Manager when you open the software. The file is divided into sections: General, Sequencing, Genotyping (STR genotyping), and SNP (SNP genotyping). Figure 11-1 shows the MegaBACE.ini file with the SNP fields highlighted.

Table 11-1 describes the fields in the SNP section of the MegaBACE.ini file. See chapter 2 for a description of the other fields in the MegaBACE.ini file.

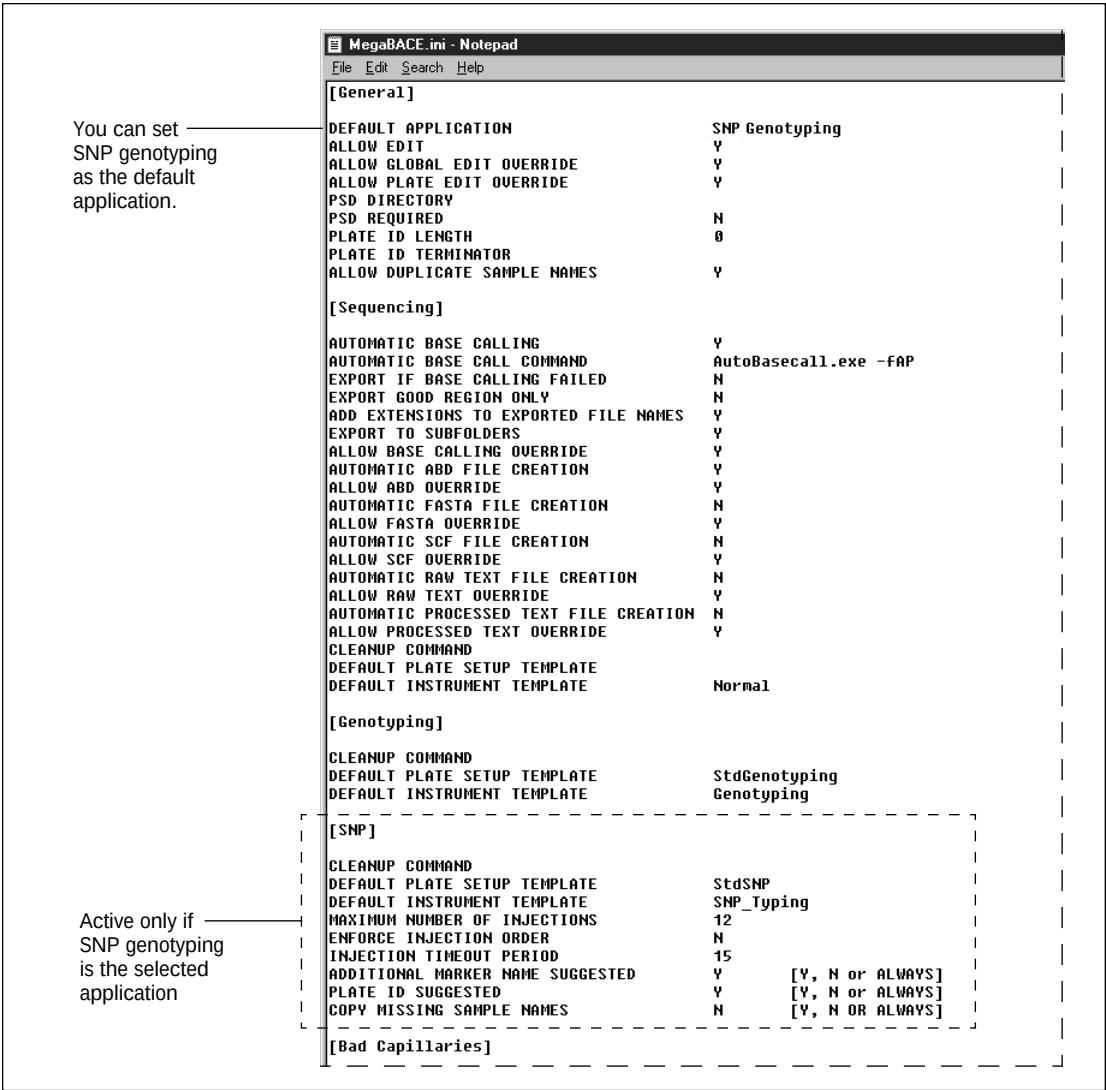


Figure 11-1. The SNP fields in the MegaBACE.ini file. Note that the file contains a section for each application, which allows you to specify a default plate setup and instrument template.

Table 11-1. The SNP genotyping fields in the MegaBACE.ini file

Field name	Description	Setting type	Default
Cleanup command	Allows the administrator to specify a program or batch file to be executed at the end of the run. This can be useful for transferring the raw data files to another workstation for analysis or archiving. See appendix C for a more detailed description.	String consisting of the path to the program or batch file (full or relative path to the instrument control folder)	none
Default plate setup template	Specifies a default plate set setup template.	String	StdSNP
Default instrument template	Specifies a default instrument control template.	String	SNP_Typing
Maximum number of injections	Specifies an upper limit and default value for the NUMBER OF INJECTIONS field that is specified in the plate set setup (section 13.2.3) for a run. If a number is not specified, the software does not impose a limit.	Integer	12
Enforce injection order	Determines whether the operator must inject the plates in the order specified in the master .psd file. Y = If the plate set .psd file specifies a PLATE SET (section 13.4.2), the bar code the operator enters in every injection must match a specific item and the order specified in the PLATE SET field. N = If the plate set .psd file specifies a PLATE SET, the bar code the operator enters in every injection must match one of the items in the PLATE SET field. The order of injections does not matter.	Y/N	N
Injection timeout period	Determines the time (minutes) the software allows for completing another injection after the injection INTERVAL TIME, specified in the plate set setup, has expired. The software displays a timer in the Additional Injection window. If the time expires before the injection is complete, the software aborts the protocol. If the Sleep After This Run protocol was selected for the run, the software also puts the instrument in sleep mode.	Integer	15

Table 11-1. The SNP genotyping fields in the MegaBACE.ini file (continued)

Field name	Description	Setting type	Default
Additional marker name suggested	Determines if a window appears when an injection is missing a SNP marker name and prompts the operator for the information: Y = The window appears if an injection is missing a SNP marker name. The operator can provide a SNP marker name. The operator can also choose to suppress the reappearance of the window and use the latest provided SNP marker name for any injections missing this attribute in the run. N = The window does not appear. For any injections missing this attribute, the software automatically copies the latest SNP marker name provided in a previous injection. ALWAYS = The window appears for every injection for which there is a missing SNP marker name. The operator cannot suppress the reappearance of the window.	Y, N, or ALWAYS	Y
Plate ID suggested	Determines if a window appears when an injection is missing a plate ID and prompts the operator for the information: Y = The window appears if an injection is missing a plate ID, but the operator can suppress the reappearance of the window for the remaining injections in the run. N = The window does not appear. ALWAYS = The window appears for every injection for which there is a missing plate ID. The operator cannot suppress the reappearance of the window.	Y, N, or ALWAYS	Y

Table 11-1. The SNP genotyping fields in the MegaBACE.ini file (continued)

Field name	Description	Setting type	Default
Copy Missing Sample Names	For an injection missing all sample names, determines if the software proceeds without sample names, prompts the operator for input, or copies the sample names from a previous injection: N = The software does not prompt the operator for input and proceeds with the injection without sample names. Y = The software prompts the operator to choose whether to copy the sample names from the latest plate containing at least one sample name or proceed with the injection without sample names. ALWAYS = The software does not prompt the operator, but automatically copies the sample names from the latest plate containing at least one sample name.	Y, N, or ALWAYS	N

11.2 Setting SNP genotyping as the default application

The operator can change the selected application in the Instrument Control Manager by choosing the Applications command from the View menu. However, when you restart the Instrument Control Manager, the software automatically resets the selected application to the default specified in the MegaBACE.ini file.

To set SNP genotyping as the default application, you use the following field and value in the MegaBACE.ini file:

Field in MegaBACE.ini file	Value
DEFAULT APPLICATION	SNP Genotyping



11.3 Setting default parameter templates for SNP genotyping

The Instrument Control Manager allows you to use templates to define the setup parameters for the plate set and the instrument control parameters for multiple runs. You can specify default templates in the MegaBACE.ini file.

If you enable the edit mode in the MegaBACE.ini file, the operator can modify the values for a given run or choose a different template in the Instrument Control Manager. However, when the Instrument Control Manager is restarted, the software automatically resets the selected template to the default specified in the MegaBACE.ini file.

11.3.1 Setting the default plate set setup template for SNP genotyping

If you specify a default template in the MegaBACE.ini file, the values from the plate set setup template appear automatically in the Plate Set Setup window. As figure 11-1 and table 11-1 show, the software has been configured to use the StdSNP.tpl template as the default for SNP genotyping. If you need to specify a different default template for the SNP genotyping application, you can use the following field in the SNP section of the MegaBACE.ini file:

Field in MegaBACE.ini file	Value
DEFAULT PLATE SETUP TEMPLATE	Type the file name of the .tpl file, without the .tpl file extension.

If the edit mode is enabled, the operator can select a different default plate set setup template from the Templates menu in the Plate Set Setup window.

Note: If the operator uses a plate setup data file (.psd) to define the plate set for a run, the values in the .psd file overwrite any duplicate parameters specified in a .tpl file.



11.3.2 Setting the default instrument control template for SNP genotyping

If you specify a default template in the MegaBACE.ini file, the values from the instrument control template appear automatically in the Instrument Control window when the operator performs a run. As figure 11-1 and table 11-1 show, the software has been configured to use the SNP_Typing.icp template as the default for SNP genotyping. If you need to specify a different default template for the SNP genotyping application, you can use the following field in the SNP section of the MegaBACE.ini file:

Field in MegaBACE.ini file	Value
DEFAULT INSTRUMENT TEMPLATE	Type the file name of the .icp file, without the .icp file extension.

If the edit mode is enabled, the operator can select a different instrument template from the Templates menu in the Instrument Control window.

11.4 Specifying limits for the injections

11.4.1 Specifying a maximum number of injections

Your laboratory might want to specify a maximum number of injections to make sure an operator does not exceed the number of injections possible for the given run conditions. All the samples in a multi-injection run must be injected into the capillaries before the samples from the first injection reach the capillary detection windows. The injection interval time (plate set definition) affects how many injections an operator can complete before the first injection reaches the capillary detection window. The time required for an injection interval is related to the length of the SNUPe primer. See the *MegaBACE 4000 Instrument User's Guide* for guidelines on designing a multi-injection run.

To specify a maximum number of injections, you use the following field in the MegaBACE.ini file:

Field in MegaBACE.ini file	Value
MAXIMUM NUMBER OF INJECTIONS	Integer (12 or less) Note: If the NUMBER OF INJECTIONS value is greater than the maximum specified by this field, the software displays an error message when the operator provides the plate set definition (section 12.2).

The software checks to see if a NUMBER OF INJECTIONS value is provided with the plate set setup parameters (tables 12-1 and 13-2). If the MAXIMUM NUMBER OF INJECTIONS value is—

- Less than the NUMBER OF INJECTIONS value, the software displays an error message.
- Greater than the NUMBER OF INJECTIONS value, the software uses the NUMBER OF INJECTIONS value as the number of injections proposed to the operator during the run.

If no NUMBER OF INJECTIONS value is provided, the software uses the MAXIMUM NUMBER OF INJECTIONS value as the number of injections proposed to the operator during the run.

11.4.2 Specifying the injection timeout period

Typically, an operator should inject the plate within 1–2 minutes after the Sample Injection window appears (figure 11-2) during the Inject Samples and Run protocol. However, you should allow sufficient time in the injection period in case an operator needs to leave the instrument temporarily during one of the injections.

To specify a time limit for the injection, you use the INJECTION TIMEOUT PERIOD field in the MegaBACE.ini file. The value you specify for the INJECTION TIMEOUT PERIOD (15 minutes, default) is the maximum time the operator has to inject a plate for any single injection in the multi-injection run.

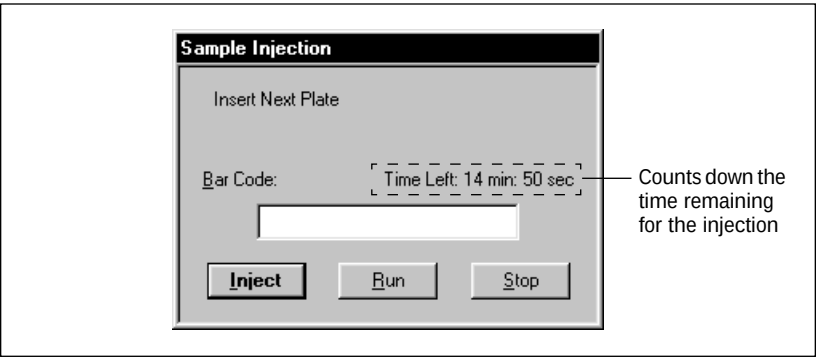


Figure 11-2. The Sample Injection window showing the timer that counts down the time remaining for the injection.

If the operator does not inject the plate before the time in the Sample Injection window elapses, the software aborts the Inject Samples and Run protocol. In addition, if the Sleep After This Run option was selected in the Instrument Control window, the software turns off the laser and the heater in the electrophoresis compartment.

11.5 Enforcing the order in which an operator injects the plates

You can use the ENFORCE INJECTION ORDER field in the MegaBACE.ini file to specify the order in which the operator injects the plates in a plate set. The value you specify for the ENFORCE INJECTION ORDER field is only applicable if you use a master .psd file to specify the list of plates in the PLATE SET field (section 13.2.2). In this case, the following values apply for the ENFORCE INJECTION ORDER field—

- **Y**—The operator must provide the plate bar codes (.psd file names) for each injection in the order specified in the PLATE SET field of the .psd file.
- **N**—The bar code the operator enters for each injection must match one of the items in the PLATE SET field, but the order of the injections does not matter.

In both cases, the software matches the bar code that the operator provides with a given item in the PLATE SET field only once. **Note:** If the NUMBER OF INJECTIONS value is greater than the number of items in the PLATE SET field (section 13.2), the operator can enter any bar codes for the remaining injections.

11.6 Prompting the operator for plate-specific information

If the plate ID, SNP marker name, or sample names are missing for an injection in a multi-injection run, you can configure the software to display a window that allows manual entry of the missing attributes (figure 11-3).

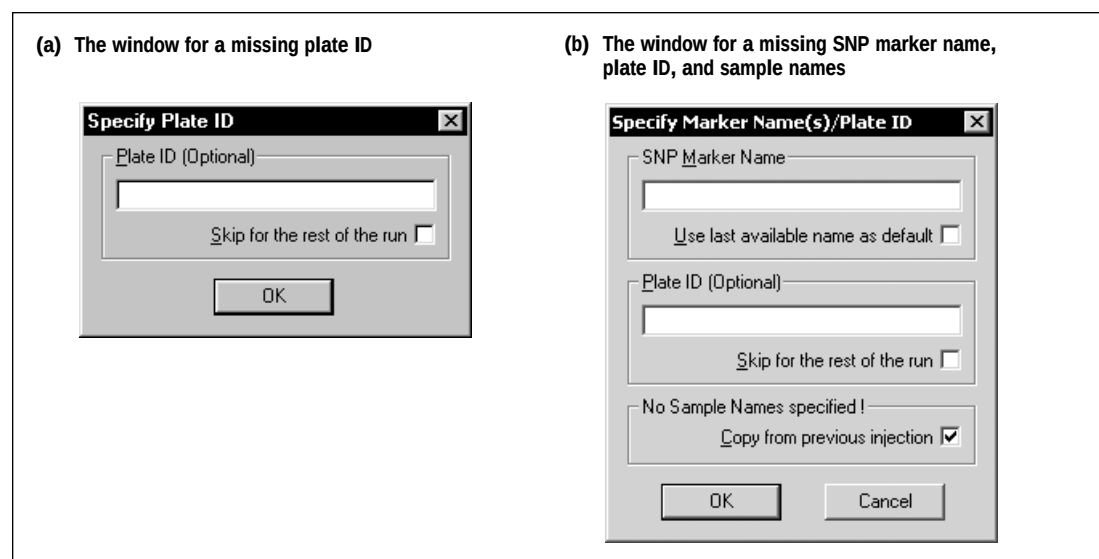


Figure 11-3. The Specify Marker Name(s)/Plate ID window changes dynamically depending on which information is missing: (a) an example of the window requesting a missing plate ID, and (b) an example of the window requesting the SNP marker name, the plate ID, and sample names.

11.6.1 Providing missing SNP marker names

The software requires at least one SNP marker name per injection. To specify one plate-level SNP marker name per injection, you can allow the operator to import the information from a .psd file for each injection or type the name in the Specify Marker Name(s) window (figure 11-3b). To specify a SNP marker name for a given well in an injection that is different than the plate-level (global) SNP marker name, the software requires a .psd file for the injection (section 13.2).

If the SNP marker name is missing for a given injection, the choices available to the operator vary depending on the value you specify in the ADDITIONAL MARKER NAME SUGGESTED field. The field is effective only if an injection is missing the plate-level marker name and does not have well-specific SNP marker names defined for all the wells in the plate. The value you specify in the ADDITIONAL MARKER NAME SUGGESTED field can depend on whether your laboratory uses—

- Different SNP markers in the same run
- The same SNP marker for multiple runs

Using different SNP marker names in the same run

If your laboratory uses multiple SNP markers in any given run, you might want the software to prompt the operator for any injections that are missing the SNP marker name. To do this, you can specify Y (default) or ALWAYS as the value for the ADDITIONAL MARKER NAME SUGGESTED field. The value that you specify depends on whether you want to allow the operator to automatically copy the last plate-level SNP marker name provided in the run for an injection that is missing the attribute.

- **Y**—Allows the operator to provide the SNP marker name and choose to automatically copy the latest plate-level SNP marker name for an injection that is missing the attribute. If an injection is missing the plate-level SNP marker name, a Specify Marker Name(s) window appears (figure 11-3b) and allows the operator to type the SNP marker name.

In addition, the operator can select the **Use last available name as default** check box. If the operator selects the check box, the Specify Marker Name(s) window does not appear again during the run. Instead, the software automatically uses the SNP marker name the operator provided for the following consecutive injections that are missing the attribute.

If the software encounters a plate in the subsequent injections that contains a plate-level SNP marker name, this latest SNP marker name is used for the subsequent consecutive injections that are missing the attribute until the software encounters an injection containing the attribute.

- **ALWAYS**—Allows the operator to provide the SNP marker name, but prevents the operator from choosing to copy the latest SNP marker name for any subsequent injections that are missing the SNP marker name.

If a Specify Marker Name(s) window appears, the **Use last available name as default** check box is unavailable. The operator must provide a SNP marker name and cannot suppress the reappearance of the window.

Using the same SNP marker name for multiple injections

If your laboratory uses the same SNP marker on many individuals, you might want to copy the SNP marker name for any injections that are missing this attribute. The default value (Y) for the ADDITIONAL MARKER NAME SUGGESTED field allows the operator to choose to copy the missing SNP marker names as described above.

Alternatively, you might want the software to copy the missing SNP marker name from the previous injection automatically instead of prompting the operator. To do this, specify N for the value of the ADDITIONAL MARKER NAME SUGGESTED field. If the first injection is missing the SNP marker name, a Specify Marker Name(s) window appears, but the window does not

reappear for subsequent injections in the same run. Instead, the software automatically uses the last plate-level SNP marker name provided during the run for the following consecutive injections that are missing the attribute. If the software encounters a plate in the subsequent injections that contains a plate-level SNP marker name, this latest SNP marker name is used for the subsequent consecutive injections that are missing the attribute until the software encounters an injection containing the attribute.

11.6.2 Providing missing plate IDs

If an injection is missing the plate ID, you can configure the software to prompt the operator for the missing plate ID or ignore the missing plate ID.

Prompting the operator for missing plate IDs

You use the PLATE ID SUGGESTED field to configure the software to display a Specify Plate ID window (figure 11-3a), which allows the operator to provide the missing plate ID. In addition, you can also specify whether the operator—

- Can suppress the reappearance of the Specify Plate ID window for the remaining injections. To enable this feature, specify Y as the value for the PLATE ID SUGGESTED field (default).

If a Specify Plate ID window appears, the operator can select the **Skip for the rest of the run** check box. If the operator selects the check box, the software proceeds without a plate ID for any injections missing this attribute.

- Cannot suppress the reappearance of the Specify Plate ID window for the remaining injections. To do this, specify ALWAYS as the value for the PLATE ID SUGGESTED field.

If a Specify Plate ID window appears, the **Skip for the rest of the run** check box is not available. A Specify Plate ID window appears for any injections in the run that are missing a plate ID.

Ignoring missing plate IDs

If you want the software to ignore any missing plate IDs and automatically proceed without prompting the operator, specify N as the value for the PLATE ID SUGGESTED field.

11.6.3 Copying sample names from a previous injection

The ability to copy all the sample names from a previous injection can be useful if your laboratory uses the same sample plate layout for multiple injections. To specify the sample names for a given injection, the software requires a .psd file (section 13.2.3). However, if all the sample names are missing for an injection,

you can allow the operator to copy the sample names from a previous injection in the same run.

You can use the COPY MISSING SAMPLE NAMES field to specify how the software handles the injections for which no sample names are specified.

Ignoring the missing sample names

If you want the software to ignore the missing sample names and proceed without prompting the operator for information, specify N as the value for the COPY MISSING SAMPLE NAMES field (default). The software uses the well ID in place of the missing sample names for the given injection.

Allowing the operator to copy the missing sample names

If you want the software to prompt the operator to copy all the sample names from the latest injection that contained sample names, specify Y as the value for the COPY MISSING SAMPLE NAMES field.

If no sample names are provided for a given injection, a window appears (figure 11-3b). In the No Sample Names Specified area, the operator can select the **Copy from previous injection** check box. If the operator selects the check box, the software copies the sample names from the latest injection containing any sample names.

Automatically copying the missing sample names

If you want the software to automatically copy all the sample names from the latest injection that contained sample names, specify ALWAYS as the value for the COPY MISSING SAMPLE NAMES field. In this case, the software does not prompt the operator to copy the missing sample names.

Chapter 12 Creating templates for SNP genotyping plate set setup parameters

For each SNP genotyping run, the Instrument Control Manager requires a plate set definition, at least one SNP marker name for each injection, and the instrument control parameters. This chapter describes plate set definitions and how to create the plate set setup templates for SNP genotyping. The topics are—

- About creating the plate set definitions (section 12.1)
- About the plate set setup parameters for SNP genotyping (section 12.2)
- Using the Plate Set Setup window to create a template (section 12.3)
- Creating a plate set setup template in a text file (.tpl) (section 12.4)
- Using a template to suppress .rsd file creation for empty capillaries (section 12.5)
- Setting and clearing the default plate set setup template (section 12.6)

For a description of the instrument control templates, see chapter 4.

12.1 About creating the plate set definitions

A plate set definition consists of the—

- **Plate set ID**—The software uses the plate set ID to name the run folder. If the plate set ID is missing, the software uses the plate ID of the first plate injected. If the plate ID is also missing, the software uses the file name of the master .psd to name the run folder.
- **Plate set setup parameters**—The run parameters that are common to all the plates in the plate set.
- **List of plates (optional)**—The list of bar codes representing the plates in the plate set. The list of plates can only be specified in a master .psd file (section 13.2.2).

12.1.1 Methods for creating the plate set definitions

The Instrument Control Manager allows you to create the plate set definitions in various ways. You can use—

- A master plate setup data file (.psd) that an operator can import at the start of the Inject Samples and Run protocol. The .psd file can specify the plate set ID, the plate set setup parameters, and the set of plates included in the multi-injection run. In addition, the master .psd file can include the plate-specific attributes for the first plate to be injected. For a description of the .psd files for SNP genotyping, see chapter 13.
- The Plate Set Setup window (figure 12-1) to create the definitions before the run. The operator can select the precreated plate set definitions at the start of the Inject Samples and Run protocol.

12.1.2 The methods for creating plate set definitions in the Plate Set Setup window

You or an operator can create the plate set definition before the run using the Plate Set Setup window (figure 12-1). To do this, you or the operator can—

- **Import a master .psd file**—You can use a master .psd file to enter the plate set definition for a single plate set. When the .psd file name is entered in the Plate Set ID box of the Plate Set Setup window, the software imports the information from the .psd file and automatically creates the plate set definition. See chapter 13 for a description of how to create .psd files.
- **Select a plate set setup template (.tpl)**—You can use a template to enter the plate set setup parameters that are visible in the Plate Set Setup window for multiple plate sets (section 12.2). In the MegaBACE.ini file, you can specify a default template that appears automatically when an operator creates a new plate set definition in the Plate Set Setup window. The software includes a default template for the SNP genotyping application (StdSNP.tpl).

A plate set setup template cannot include fields for the plate set ID, the SNP marker name, or the list of plates in the run. However, you can use a master .psd file in addition to the template to specify these additional parameters.
- **Manually enter or modify the setup parameters**—To manually enter or change parameters in the Plate Set Setup window, you must enable the edit mode in the MegaBACE.ini file. See section 2.3 for a description of the edit mode.

See the *MegaBACE 4000 Instrument User's Guide* for a description of how to create plate set definitions in the Plate Set Setup window.

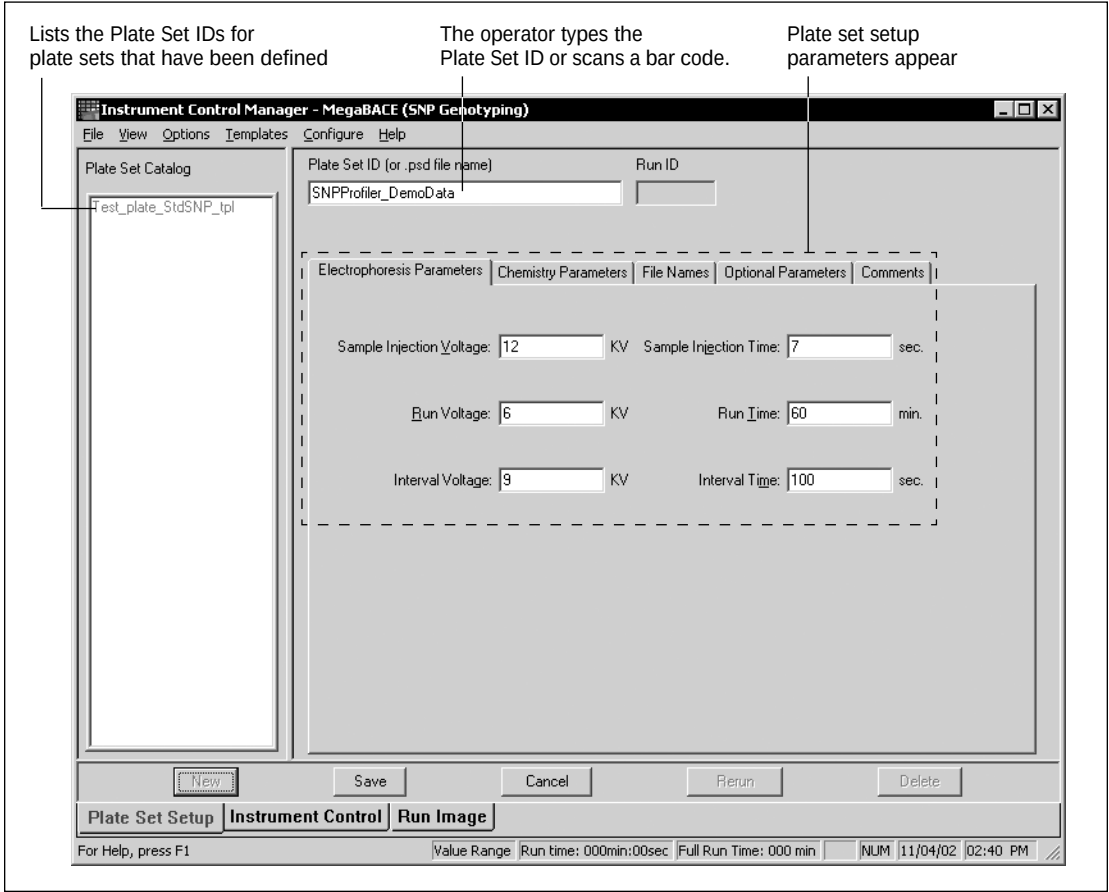


Figure 12-1. The Plate Set Setup window displays the parameters shared by all the plates in the plate set.

12.2 About the plate set setup parameters for SNP genotyping

For a SNP genotyping run, the plate set setup parameters are common to the entire set of plates injected during the multi-injection run. The plate set setup parameters consist of the electrophoresis parameters, the chemistry parameters, the file names, plate set comments, and other parameters shared by all the plates injected in the run. The values for these parameters can be imported from a .psd file, selected from a plate set setup template, or manually entered in the Plate Set Setup window.

Table 12-1 lists the plate set setup parameters. The software requires values for all the parameters, except the optional parameters and comments.

Table 12-1. The plate set setup parameters

Field name	Value	Default value for SNP genotyping
Electrophoresis Parameters tab		
Sample injection voltage	Integer or fraction (kV). The range is 1–15 kV.	12 kV
Sample injection time	Integer (s). The range is 1–600 s.	7 s
Run voltage	Integer or fraction (kV). The range is 1–15 kV.	6 kV
Run time	Integer (min). The range is 1–720 min.	60 min
Interval voltage	Integer or fraction (kV). The range is 1–15 kV.	9 kV
Interval time *	Integer (s). The range is 1–1800 s.	100 s
Chemistry Parameters tab		
Chemistry	String. The name of a chemistry parameter set. Note: The name must match a name in the Chemistry.ini file (chapter 3) for the SNP genotyping application. Otherwise, the software will not recognize the name.	SNuPe Terminators
File Names tab		
Well_ID	String. The name of the raw sample data file (.rsd) the software produces for each well in a plate. The software uses the well ID as the file name unless alternative file names are provided in the master .psd file, a plate set setup parameters template, or the Plate Set Setup window. Note: Unlike sequencing and STR genotyping, to specify the sample names for a SNP genotyping run, you must use .psd files.	Well IDs

Table 12-1. The plate set setup parameters (continued)

Field name	Value	Default value for SNP genotyping
Optional Parameters tab		
Note: The Run temperature and PMT voltages are typically defined as part of the instrument control parameters (chapter 4). The values can be defined here for a special run. If defined here, the values temporarily overwrite the same values defined in an .icp file or the Instrument Control window.		
Run temperature	Integer (°C). The range is 27–44 °C.	No default specified
PMT1 voltage	Integer or fraction (kV). The range is 500–900 V.	No default specified
PMT2 voltage	Integer or fraction (kV). The range is 500–900 V.	No default specified
PMT3 voltage	Integer or fraction (kV). The range is 500–900 V.	No default specified
PMT4 voltage	Integer or fraction (kV). The range is 500–900 V.	No default specified
Number of injections*	Integer. The range is 1–12. The MegaBACE.ini file specifies an upper limit for this field (section 11.4.1). Note: This value specified the number of injections the software prompts the operator for during a multi-injection run.	No default specified
Comments tab		
Comments	String, including spaces, representing the comments that you want to include about the plate set.	No default specified

* See the *MegaBACE 4000 Instrument User's Guide* for guidelines on determining the interval time and the number of injections.

12.3 Using the Plate Set Setup window to create a template

You can use the Plate Set Setup window (figure 12-1) to create a plate set setup template. To do this—

1. In the Plate Set Setup window, click **New**.
2. In the following tabs, enter the parameters that you want to use for the template. See table 12-1 for a description of the parameters.
 - **Electrophoresis Parameters tab**—Type the values for the parameters in each of the boxes.
 - **Chemistry Parameters tab**—Select the name of the chemistry parameters from the **Chemistry Name** list.

- **File Names tab**—By default, the software uses the well IDs to name the files (for example, A01.rsd). If you want to use file names other than the well IDs, you can type the names in the File Names tab.

If you want to use the File Names tab to specify empty wells for which you want to suppress the raw sample data file (.rsd) creation, see section 12.5.
- **Optional Parameters tab**—(Optional) If you want to overwrite the Run Temperature and PMT Voltage values, which are provided in the instrument control parameters, you can type the values in this tab. You can also type a value for the Number of Injections.
- **Comments tab**—You can type any comments about the plate set.

See the *MegaBACE 4000 Instrument User's Guide* for details on using the Plate Set Setup window.

3. On the **Templates** menu, point to **Plate Setup Templates**, and click **Save Template**. The Save Template window appears.
4. Type the template name, and click **Save**. The template is saved in the ...\\MegaBACE\\Instrument Control\\Templates folder. The software adds the .tpl file extension.

12.4 Creating a plate set setup template in a text file (.tpl)

Important

When you create a plate setup template, you must type the .tpl file extension and save it in the ...\\MegaBACE\\Instrument Control\\Templates folder. Otherwise, the Instrument Control Manager will not find the file.

To create a plate set setup template—

1. Open a new document or an existing .tpl file using a text editor (such as Notepad) or Excel.
2. Type the field names in the first column (field column), and type the values for each field in the second column (value column). Make sure the field column and value column are separated by a tab. See table 12-1 for a list of the recognized fields for the plate set setup parameters. Figure 12-2 shows an example of the StdSNP.tpl file for SNP genotyping, which is included with the software.
3. Save the file in an appropriate format. If you use Excel, save the file as a tab-delimited file. You must type the .tpl file extension.

Note: If you enter the optional parameters listed in table 12-1, the values will overwrite the values for the same parameters defined in the Instrument Control window.

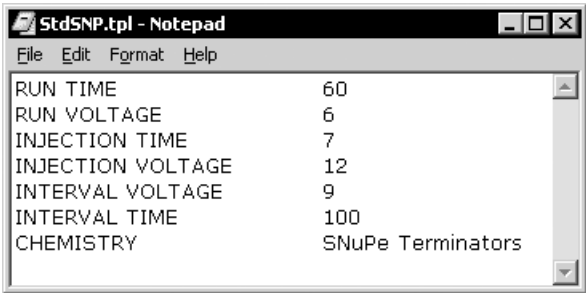


Figure 12-2. An example of a plate set setup template file (.tpl) for SNP genotyping.

Note: To specify the list of plates in a plate set before the run, the software requires the PLATE SET field in a plate set .psd file (section 13.2.2).

12.5 Using a template to suppress .rsd file creation for empty capillaries

To specify an empty capillary, you designate the corresponding well as unused. The software does not create a raw sample data file (.rsd) for any wells that are specified as not used. In the Run Image window (figure 2-2), an unused well is indicated by a yellow arrow that appears above the capillary button.

Note: To specify bad capillaries, see section 2.6.

You can specify unused wells in either a—

- **Plate setup data file (.psd)**—You can use a .psd file to specify the unused wells if the unused wells occur for a single plate (chapter 13).
- **Plate setup template**—You can use a plate setup template if the unused wells occur for multiple plates.

You can specify unused wells in a plate setup template by typing **NOT USED** (case insensitive) as the sample name specified in either the File Names tab of the Plate Set Setup window or the .tpl file (section 12.4).

To specify an unused well in the File Names tab of the Plate Set Setup window—

1. In the **File Names** tab of the Plate Set Setup window (figure 12-3), adjust the width of the column that includes the unused well. To increase the column width, click the **RIGHT ARROW** key in the File Names tab. To decrease the column width, click the **LEFT ARROW** key in the File Names tab.
2. Type **NOT USED** (case-insensitive) as the sample name for the unused well.

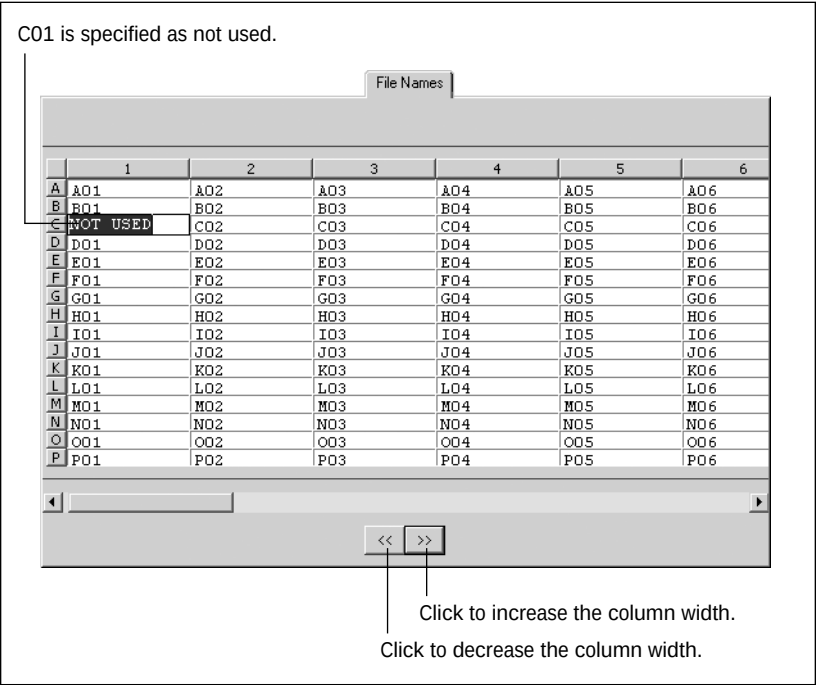


Figure 12-3. The File Names tab of the Plate Set Setup window.

12.6 Setting and clearing the default plate set setup template

You can set a default plate set setup template in the—

- **SNP section of the MegaBACE.ini file**—To set a default plate set setup template for SNP genotyping in the MegaBACE.ini file, see section 11.3.1.
- **Templates menu in the Plate Set Setup window**—The default template is active only until you close the Instrument Control Manager. When the software is restarted, the default plate set setup template specified in the SNP section of the MegaBACE.ini file is active until you change it again.

12.6.1 Using the Templates menu to set a default template

To set the default plate set setup template—

1. In the Instrument Control Manager, click the **Plate Set Setup** tab to display the Plate Set Setup window (figure 12-1).
2. On the **Templates** menu, point to **Plate Setup Template** and click **Set Default**. The Set Default window appears.
3. Select the name of the template you want to use as the default, and click **Open**. The Instrument Control Manager designates the template you selected as the default template.

The next time you click New in the Plate Set Setup window, the values from the default template appear in the parameter boxes.

12.6.2 Clearing the default plate set setup template

To clear the default plate set setup template—

1. In the Instrument Control Manager, click the **Plate Set Setup** tab to display the Plate Set Setup window (figure 12-1).
2. On the **Templates** menu, point to **Plate Setup Templates**, and click **Clear Default**. The template you cleared is no longer the default template.

When you click New in the Plate Set Setup window, no values appear in the parameter boxes. You can select another default template (section 12.6.1), or use the Plate Set Setup window to create another template, save it, and then select it as the default template.

Important

When you restart the software, the default templates specified in the MegaBACE.ini file are active.

Chapter 13 Creating plate setup data files (.psd) for SNP genotyping

For each SNP genotyping run, the Instrument Control Manager requires a plate set definition, at least one SNP marker name for each injection, and the instrument control parameters. Although not required, plate setup data files (.psd) are recommended to automate how an operator provides the plate set definition and plate-specific attributes during a multi-injection run. This chapter describes how to create .psd files and provides examples of the plate set setup templates and instrument parameter templates for SNP genotyping. The topics are—

- Creating plate setup data files (.psd) for SNP genotyping (section 13.1)
- Specifying the plate set definition in a master .psd file (section 13.2)
- Using a .psd file to suppress .rsd file creation for empty capillaries (section 13.3)
- Specifying the attributes for a plate in a .psd file (section 13.4)

13.1 Creating plate setup data files (.psd) for SNP genotyping

13.1.1 About the types of .psd files

For SNP genotyping you can use two types of .psd files, one as the master file for the plate set and the first plate you inject and one for each subsequent plate in the plate set—

- **Master .psd file**—Figure 13-1 shows a master .psd file. The master .psd file for the plate set can list the set of plates included in a multiple-injection run and define the plate setup parameters shared by all the plates in the run. The master .psd file can also include the plate-attribute fields for the first plate that you inject. For example, you can specify a plate ID and the SNP marker name for the first plate, and you can associate a sample name with each well in the first plate.

Note: If you want to define the sample names for a plate in a multi-injection run, a .psd file is required.

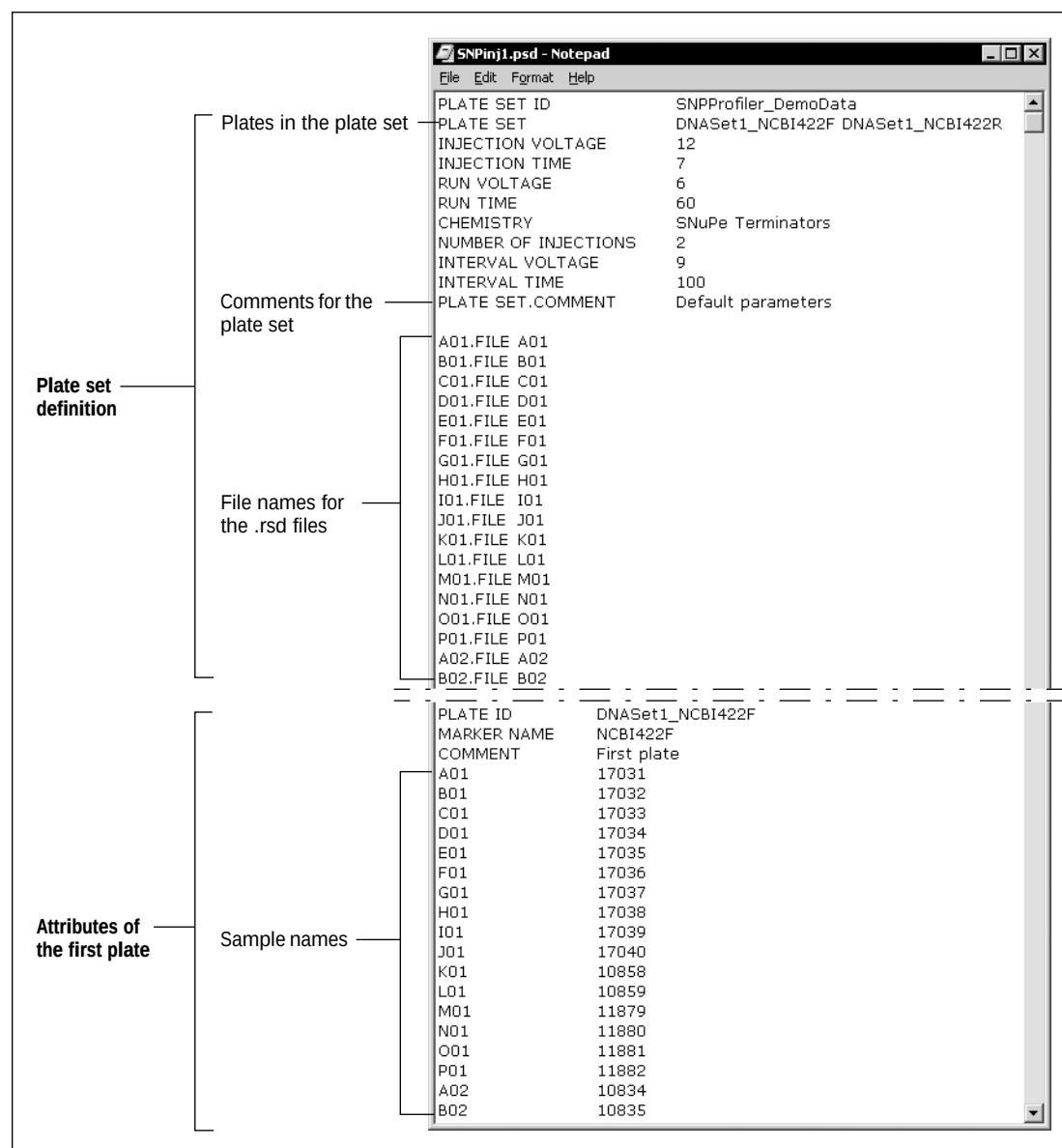


Figure 13-1. Example of master .psd file that defines both the plate set and the first plate to be injected.

- **Plate-only .psd files**—For each plate in the plate set, you can use a plate-only .psd file to specify the plate-level and well-level information for a given plate (figure 13-2). If you want to associate a sample name with each well in a plate, you must use a .psd file to define the sample names.

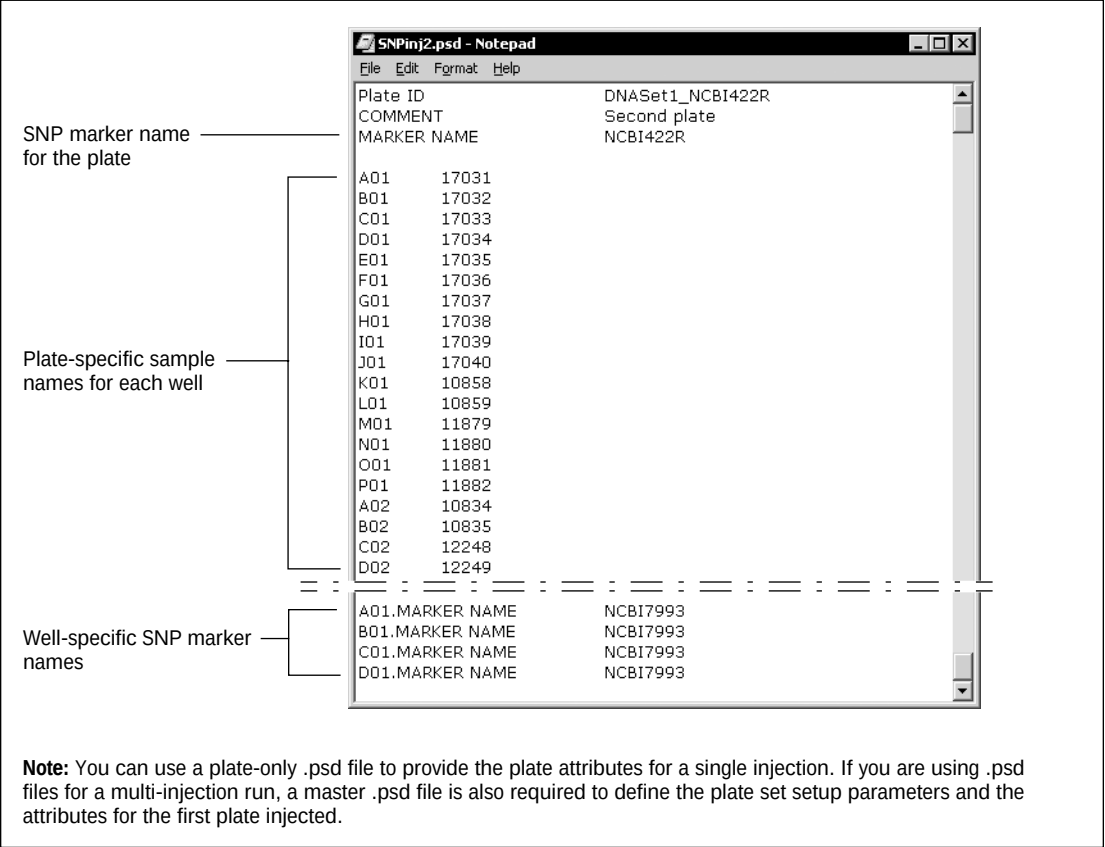


Figure 13-2. Example of a plate-only .psd file defining one plate (injection) in the plate set.

13.1.2 Creating a plate setup data file (.psd)

The .psd files can be created using any software that supports tab-separated text. When the .psd file is imported, the Instrument Control Manager stores the information in the extended header of the raw sample data files (.rsd) for all the files in the run. The well-level information, such as sample name, becomes part of the extended header in the .rsd file for the well.

Note: All fields in a .psd file are optional. However, if the required plate set setup parameters are not provided in a .psd file, the Plate Set Setup window must be used to create the plate set definition (section 12.1).

To create a .psd file—

1. Open a new document or an existing .psd file using software such as Notepad or Excel. **Note:** The MegaBACE software includes examples of .psd files that you can use as templates.
2. Using a two-column format, type the appropriate fields and values for the .psd file. The fields and values that you can use depend on whether you are creating a—

- **Master .psd file**—Figure 13-1 shows an example of a master .psd file that includes the plate set definition. Table 13-1 lists the fields for the master .psd file. In addition to these fields, the master .psd file can contain fields that define the attributes of the first plate injected in the run (table 13-2).

Note: If both a .psd file and template are used to define the same plate set setup parameter, the values provided in the file imported last overwrite the values for any identical fields.

- **Plate-only .psd file**—Figure 13-2 shows an example of a .psd file that defines a single plate (injection). Table 13-2 lists the plate-attribute fields you can use in a .psd file for a single plate.

Make sure you type the field names in the first column (the field column), and the values for each field name in the second column (the value column). Separate the field column from the value column using a tab.

3. Save the file in an appropriate format—
 - **File name**—Make sure the file name is an exact match to the plate set ID or bar code that the operator will use to import the file. Make sure you use the .psd file extension.

- **Storage location**—Make sure you store the .psd file in the ...\\MegaBACE\\Psd folder (or the folder you specified as the Psd directory in the MegaBACE.ini file); otherwise, the Instrument Control Manager cannot find the file.

Note: If you use Excel, save the file as a tab-delimited file.

Table 13-1. The fields in the master .psd file for a plate set

Field name	Value
Plate set ID	String consisting of any characters except spaces or other characters not allowed in file names.
Plate set	(Optional) String consisting of bar codes separated by spaces. The bar codes must match the file names of the .psd files in the plate set.
User-defined plate set fields	Fields that can be imported, but are not used by the software. The software stores the information from these fields in the extended header of the files.
Fields that can be specified in either a master .psd file or a plate set setup template (chapter 12):	
Injection voltage	Integer or fraction (kV). The range is 1–15 kV.
Injection time	Integer (s). The range is 1–600 s.
Run voltage	Integer or fraction (kV). The range is 1–15 kV.
Run time	Integer (min). The range is 1–720 min.
Interval voltage	Integer or fraction (kV). The range is 1–15 kV.
Interval time	Integer (s). The range is 1–1800 s.
Chemistry	String. The name of a chemistry parameter set. Note: The name must match a name in the Chemistry.ini file (chapter 3) for the SNP genotyping application. Otherwise, the software will not recognize the name.
Optional parameters: The run temperature and PMT voltages are typically defined as part of the instrument control parameters (chapter 4). If the values are defined here for a special run, the values temporarily overwrite the same values defined in an .icp file or the Instrument Control window.	
Run temperature	Integer (°C). The range is 27–44 °C.
PMT1 voltage	Integer (kV). The range is 500–900 V.
PMT2 voltage	Integer (kV). The range is 500–900 V.
PMT3 voltage	Integer (kV). The range is 500–900 V.
PMT4 voltage	Integer (kV). The range is 500–900 V.
Number of injections	Integer Note: The upper limit can be specified in the SNP section of the MegaBACE.ini file (section 11.4).

Table 13-1. The fields in the master .psd file for a plate set (continued)

Field name	Value
Well_ID.FILE (where well_ID is the location of the well, such as A01)	String representing the file name that you want to assign the designated well location. If you use spaces in the file names, the Instrument Control Manager converts the spaces to underscore characters in the file name. Note: An .rsd file for a given well location contains the sample data from all the plates injected during the run.
Plate Set.comment	String, including spaces, representing the comments for the plate set.

Table 13-2. The plate-attribute fields used in any .psd file for SNP genotyping *

Field name	Value
Plate ID	String representing the plate ID or bar code for the plate.
Marker name	String representing the SNP marker name for the plate.
Comment	String, including spaces, representing comments about the plate or injection.
Well_ID or well_ID.NAME (where well_ID is the location of the well, such as A01)	String representing the sample name that you want to assign to the designated well location. If you use spaces in the sample names, the Instrument Control Manager converts the spaces to underscore characters in the file name. Note: The software does not use this label as a file name because the sample name is specific to a well for a single plate in the plate set.
Well_ID.MARKER NAME (where well_ID is the location of the well, such as A01)	String representing the SNP marker name that you want to assign to the designated well location.
User-defined fields	Fields that can be imported, but are not used by the software. The software stores the information from these fields in the extended header of the files.

* To define plate attributes such as sample names, you must use a .psd file for the injection. If you use a .psd file for the first injection, the .psd file typically includes the plate set setup parameters as well as the attributes for the first plate.

13.2 Specifying the plate set definition in a master .psd file

This section describes the fields you can use in a master .psd file for a SNP genotyping run. See section 13.1 for instructions on how to create the .psd file.

13.2.1 Specifying the plate set ID


To define the plate set ID, you can use the PLATE SET ID field in the .psd file. The software uses the plate set ID to name the run folder. If the PLATE SET ID field does not contain a value, the software uses the value provided for the PLATE ID field. If the PLATE ID is also missing, the software uses the .psd file name (bar code) as the plate set ID and to name the raw run folder.

When an operator enters the text or bar code in the Plate Set ID area of the Plate Set Setup window or the Select a Plate Set window, the Instrument Control Manager looks for a .psd file with a file name that matches the text entered and automatically imports the plate set setup parameters that are contained in the .psd file.

Important

Regardless of the method used to enter the plate ID, the text entry or bar code must match the name of the .psd file in order for a .psd file to be imported.

Figure 13-3 shows an example of a master .psd file and the Plate Set Setup window.

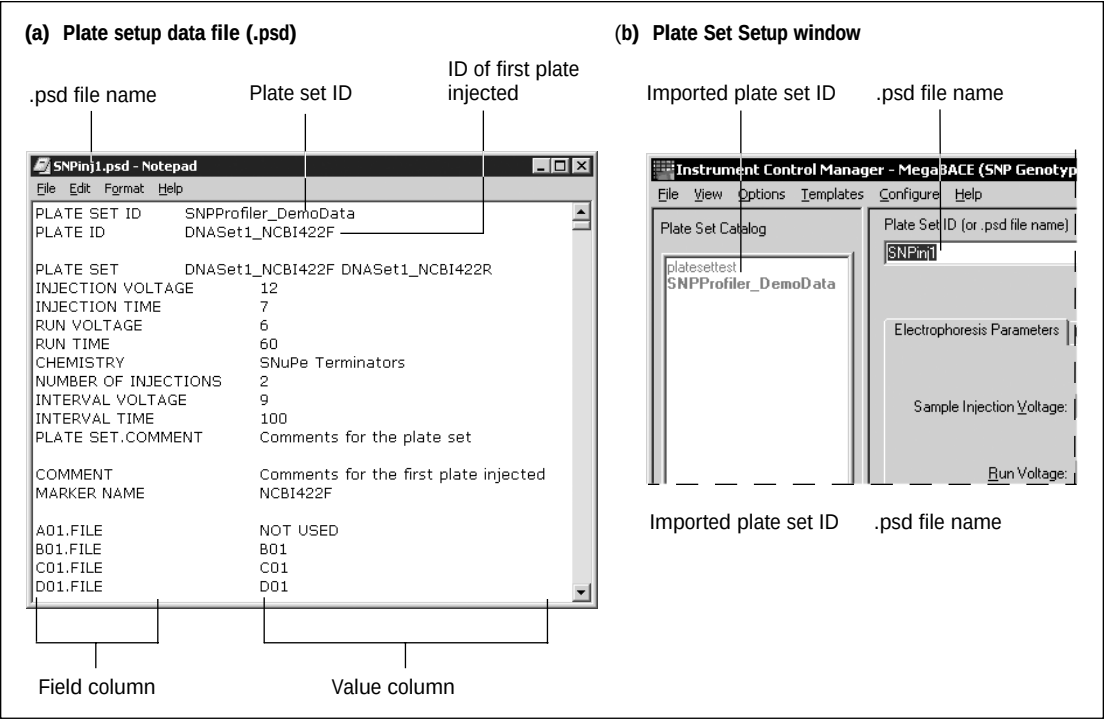

The values you use for the PLATE SET field should agree with the values you specify for the following fields:

- NUMBER OF INJECTIONS field in the .psd file. If the NUMBER OF INJECTIONS value is—
 - Less than the number of bar codes in the PLATE SET, the NUMBER OF INJECTIONS value determines the number of injections. During each injection, the operator must provide one of the bar codes in the PLATE SET.
 - Greater than the number of bar codes in the PLATE SET, the operator can use any bar codes (.psd file names) for the injections remaining after using all the bar codes listed for the PLATE SET.
- ENFORCE INJECTION ORDER field in the SNP section of the MegaBACE.ini file (section 11.5), which allows you to enforce the order in which the operator injects the plates listed in the PLATE SET.

13.2.3 Defining the number of injections

To define the number of injections, you can use the NUMBER OF INJECTIONS field in the .psd file. The value you specify for this field should agree with the values specified for the following fields:

- PLATE SET field in the .psd file. See section 13.2.2 for details.
- MAXIMUM NUMBER OF INJECTIONS field in the MegaBACE.ini file (section 11.4.1). If the NUMBER OF INJECTIONS value exceeds the maximum number specified in the MegaBACE.ini file, the software displays an error message when the operator tries to import the .psd file.

13.3 Using a .psd file to suppress .rsd file creation for empty capillaries

You can indicate unused wells by typing **NOT USED** (case-insensitive) as the file name specified in the .psd file (figure 13-3a). For every unused well, a yellow arrow appears above the capillary button in the Run Image window (figure 2-2). The software does not create an .rsd file for a well specified as not used. For details on suppressing .rsd files for defective capillaries, see section 2.6.

13.4 Specifying the attributes for a plate in a .psd file

This section describes the fields you can use in a .psd file to describe the attributes of a single plate (injection) in a SNP genotyping run. The software does not display the plate-attribute fields in the Plate Set Setup window, because the fields apply to a single plate in the plate set. You can use the plate-attribute fields in both the master .psd file and the plate-only .psd files.

Note: If you do not use a master .pdf file to define the setup parameters for the plate set and the attributes for the first plate, you cannot use the plate-only .psd files to define the plate attributes for each injection.

Figure 13-2 shows an example of a plate-only .psd file. Table 13-2 lists the fields you can use in a .psd file to define the attributes of a specific plate in the plate set. These fields can be defined only in a .psd file. See section 13.1 for instructions on how to create the .psd file.

13.4.1 Specifying the plate ID

In the .psd file for a given plate, you can specify the plate ID. Alternatively, you can configure the MegaBACE.ini file so that the software always prompts the operator for the plate ID for each injection (section 11.6).

13.4.2 Specifying the SNP marker name

The Instrument Control Manager requires at least one SNP marker name to be specified for each injection in a SNP genotyping run. Multiple SNP marker names for each plate (injection) are allowed. To specify the names of the SNP markers for a plate, you can use the following fields in the .psd file:

- **MARKER NAME**—Specifies a SNP marker name for the entire plate.
- **Well_ID.MARKER NAME**—Specifies the SNP marker name for the selected well ID (figure 13-2). The well-specific SNP marker name overwrites the plate SNP marker name for the given well.

Alternatively, if you do not use .psd files, you can configure the MegaBACE.ini file (section 11.6) so that the software always prompts the operator for the SNP marker name for each injection during the run. This method allows the operator to enter only the global SNP marker name for each plate.

Note: SNP Profiler requires a profile for each SNP marker name found in the extended header of the .rsd files. See the *MegaBACE SNP Profiler User's Guide* for a description of the SNP Marker Editor.



13.4.3 Specifying sample names for a given plate

To specify a sample name for a given well on a plate in a SNP genotyping run, you must use the well_ID.NAME field or well_ID field in a .psd file.

If your laboratory routinely uses the same sample plate layout for multiple injections, you might want to copy the sample names from one injection to the next. See section 11.6.3 for a description of how to configure this option.

Part five

Protocol modification and file storage

Chapter 14 Bypassing or overriding a protocol

This chapter describes how to bypass or override a protocol. The topics are—

- Bypassing the selected protocol (section 14.1)
- Overriding the contents of a protocol (section 14.2)

In a typical operation, the software selects the next protocol in the workflow. You might need to use a protocol that is different from the one selected in the List of Protocols in the Instrument Control window (figure 14-1).

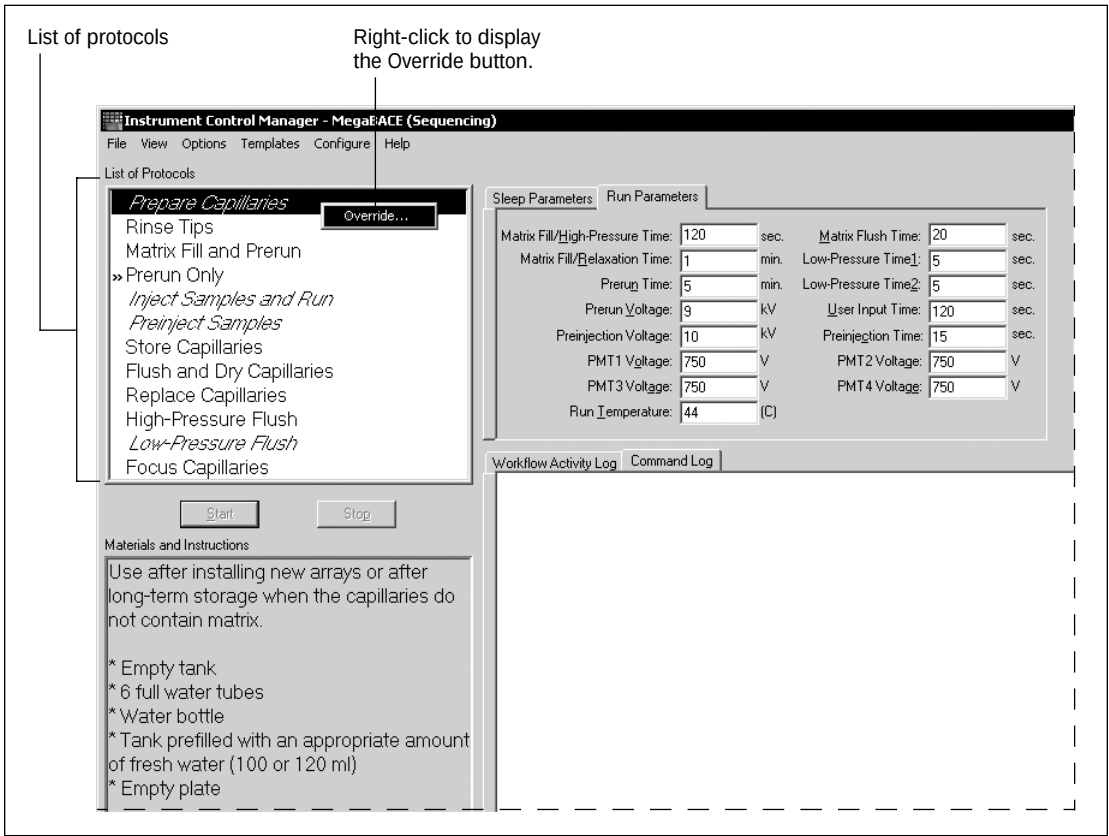


Figure 14-1. The List of Protocols in the Instrument Control window.

The Instrument Control Manager allows you to—

- Bypass the protocol the software selects. You can select any protocol, even if it is disabled (section 14.1).
- Override all or part of the protocol (section 14.2).

14.1 Bypassing the selected protocol

To bypass the selected protocol—

1. In the Instrument Control window, select the protocol you want to run, and click the **right mouse button**. The Override button appears (figure 14-1).
2. Click the **Override** button. A window appears and asks for a password (figure 14-2).



Figure 14-2. The Enter Password window.

3. Type the appropriate password and click **OK**. The Override Commands window appears (figure 14-3).

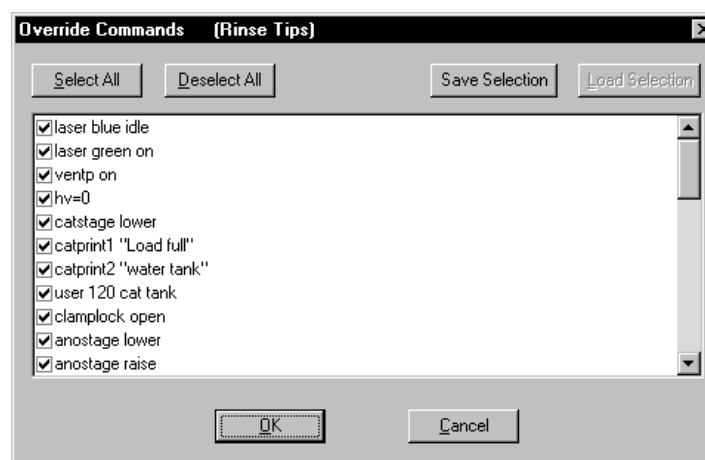


Figure 14-3. The Override Commands window.

4. Click **OK** to run all the commands for the selected protocol. The window closes. In the Instrument Control window (figure 14-1), the selected protocol name appears bold to show that it is running.

14.2 Overriding the contents of a protocol

You can override all or part of the protocol the software selects as the next logical protocol. To do this—

1. With the protocol selected, click the **right mouse button**. The Override button appears (figure 14-1).
2. Click the **Override** button. A window appears and asks for a password (figure 14-2).
3. Type the appropriate password and click **OK**. The Override Commands window appears (figure 14-3).
4. In the Override Commands window—
 - Select the commands you want to run. You can click **Deselect All** and select the desired commands. Alternatively, you can deselect only the commands you do not want to run.
 - (Optional) Click **Save Selection** to save a specific set of commands you want to re-select later. You can save one set of selected commands for each protocol. To re-select a previously saved set of selected commands, click **Load Selection**.
 - Click **OK** to close the window.

In the Instrument Control window (figure 14-1), the selected protocol name appears bold to show that it is running.

Chapter 15 File storage

This chapter describes—

- Raw sample data storage (section 15.1)
- Base-called sample data storage (sequencing only) (section 15.2)
- System initialization file storage (section 15.3)
- Parameter configuration file storage (section 15.4)

15.1 Raw sample data storage

This section describes how the sample data files are stored and what the files contain.

15.1.1 Raw sample data file content

For each run, the Instrument Control Manager creates a raw sample data file (.rsd) for the data collected from each well. Each .rsd file contains—

- Plate ID (sequencing and STR genotyping) or plate set ID (SNP genotyping)
- Well location
- If provided, sample name (sequencing and STR genotyping) or sample names (SNP genotyping)
- Run data: run ID, instrument ID, instrument control parameters, plate setup parameters, and date
- Information about the electrical current
- Raw electropherogram (signal intensities)

15.1.2 Naming and storage conventions (sequencing and STR genotyping only)

By default, the software uses the sample names that you provide as the file names for the .rsd files. If the sample names that you assign for a given plate are not unique, the Instrument Control Manager appends the well ID to the file name (sample_name_well_ID) only if the ALLOW DUPLICATE SAMPLE NAMES field in the MegaBACE.ini file is set to Y. If no sample names are provided, the software uses the well ID as the file name for the .rsd files.

Alternatively, you can use a plate setup data file (.psd) to specify the explicit file names for the .rsd files. For sequencing, see section 7.5. For STR genotyping, see section 10.4.

Each time you run a plate on the MegaBACE instrument, the Instrument Control Manager creates a raw run folder that includes the .rsd files for the run. The MegaBACE system uses the plate ID and the run ID to name the raw run folders (plateID_runID). You assign the plate ID when you set up the plate in the Instrument Control Manager software. The run ID defaults to Run01. Table 15-1 summarizes the conventions for naming and storing the raw sample data files and folders.

Table 15-1. Naming and storage conventions for the raw sample data (sequencing and STR genotyping only)

File or folder type	Name	Location
Raw sample data file (.rsd)	Sample_name.rsd, well_ID.rsd, or a combination of both	PlateID_runID (raw run folder)
Raw run folder	PlateID_runID	...\MegaBACE\Data (unless you select a different location)

The File Storage command on the Options menu allows you to change the storage location of the data at any time. This feature only affects the data storage of subsequent runs. It does not change the location of the data already stored on the computer. See the *MegaBACE 4000 Instrument User's Guide* for a description of how to change the file storage location.

15.1.3 Naming and storage conventions (SNP genotyping only)

By default, the software uses the well ID as the .rsd file name, for example, A01.rsd. Alternatively, you can use a master .psd file, a plate set setup template, or the Plate Set Setup window to specify the file names for the .rsd files (chapter 13).

Note: Unlike sequencing and STR genotyping, to specify the sample names for a SNP genotyping run, you must use .psd files.

Each time you perform a multi-injection run, the Instrument Control Manager creates a raw run folder that includes the .rsd files for the run. The software uses the plate set ID and the run ID to name the raw run folders (plate set ID_runID). If the plate set ID is missing, the software uses the plate ID of the first plate injected. If the plate ID is also missing, the software

uses the file name of the master plate setup data file (.psd) to name the run folder. The run ID defaults to Run01. Table 15-2 summarizes the conventions for naming and storing the raw sample data files and folders.

Table 15-2. Naming and storage conventions for the raw sample data (SNP genotyping only)

File or folder type	Name	Location
Raw sample data file (.rsd)	Well_ID.rsd	Plate set ID_runID (raw run folder)
Raw run folder	Plate set ID_runID	...\MegaBACE\Data (unless you select a different location)

The File Storage command on the Options menu allows you to change the storage location of the data at any time. This feature only affects the data storage of subsequent runs. It does not change the location of the data already stored on the computer. See the *MegaBACE 4000 Instrument User's Guide* for a description of how to change the file storage location.

15.2 Base-called sample data storage (sequencing only)

You can specify that the software create a subfolder for each export file format that you use after automatic base calling (chapter 5). Sections 15.2.1 and 15.2.2 describe the naming and storage conventions. See the *MegaBACE 4000 Instrument User's Guide* for a description of how to change the storage location for the base-called data.

15.2.1 Naming and storage conventions for the base-called and exported files

Table 15-3 provides an overview of the naming and storage conventions for the base-called and exported files. Note that the base-called sample data files obtain the file names from the .rsd files.

Table 15-3. Naming and storage conventions for the base-called and exported files

File type	File contents	File name	Folder storage
Base-called	Plate ID, well location, sample name, instrument ID, base caller name, called sequence, sequence starting and ending points, read length for the high-quality data, quality values, electrical current information, base-called electropherogram for the sample.	Sample_name.esd, well_ID.esd	PlateID_basecallerID (analyzed run folder)
Exported—ABD	All the raw data, including electropherogram and spectral separation table. If base-called, contains all the raw and base-called data, including the called sequence.	Sample_name.abd, well_ID.abd	ABD subfolder (default) Note: In the MegaBACE.ini file, you can specify that the files be stored at the top level of the run folder (table 5-1 and chapter 5).
Exported—SCF	Raw data and the electropherograms. If base-called, includes the raw and base-called data, and the called sequence.	Sample_name.scf, well_ID.scf Note: The .scf extension might not display in Windows Explorer.	SCF subfolder (default) Note: In the MegaBACE.ini file, you can specify that the files be stored at the top level of the run folder (table 5-1 and chapter 5).
Exported—FASTA	File name, file location, and the called sequence for each well in the base-called plate (default). You can use the MegaBACE.ini file to specify other options (section 5.5.2): • One .seq file for each .esd file • One .seq file for each base-called plate • One .seq file for all base-called plates	.seq	PlateID_basecallerID folder (default)
Exported—ASCII text	Plate name, sample name, instrument ID, base caller name, sequence starting and ending points, electrical current information, scan rate, instrument parameters, the number of bases called.	Sample_name.txt, well_ID.txt	Text subfolder (default) Note: In the MegaBACE.ini file, you can specify that the files be stored at the top level of the run folder (table 5-1 and chapter 5).

15.2.2 Naming and storage conventions for the base-called run folders

The different types of base-called data (section 15.1.2) are stored in different folders. Table 15-4 lists the run folders, exported data folders, and the naming conventions for these folders.

The software uses the plate ID, the base caller ID, the base caller version number (n_nnn), and the base caller method to name the analyzed run folders. The base caller IDs are—

- M = Molecular Dynamics™
- Cp = All Cimarron™ base callers

Table 15-4. Naming conventions for the analyzed run folders

Folder type	Folder name	Storage location
Analyzed (Molecular Dynamics base caller)	PlateID_Mn_nnn_MD1	...\MegaBACE\AnalyzedData folder (default)
Analyzed (Cimarron 1.53 Phat base caller)	PlateID_Cp1_53000_MD1	...\MegaBACE\AnalyzedData folder (default)
Analyzed (Cimarron 1.53 Slim Phredify base caller)	PlateID_Cp1_53000_sl_ph_MD1	...\MegaBACE\AnalyzedData folder (default)
Analyzed (Cimarron 3.12)	PlateID_Cp312_MD1	...\MegaBACE\AnalyzedData folder (default)
Analyzed (Cimarron 3.12 Even Spacing)	PlateID_Cp312_es_MD1	...\MegaBACE\AnalyzedData folder (default)
Analyzed (Cimarron 3.12 Aligned base caller)	PlateID_Cp312_a_MD1	...\MegaBACE\AnalyzedData folder (default)
Exported ABD	Abd	PlateID_basecallerID folder (default)
Exported SCF	Scf	PlateID_basecallerID folder (default)
Exported text	Text	PlateID_basecallerID folder (default)

15.3 System initialization file storage

Figure 15-1 shows the system initialization file and parameter configuration file locations. Table 15-5 lists the system initialization files that you can modify to configure the Instrument Control Manager.

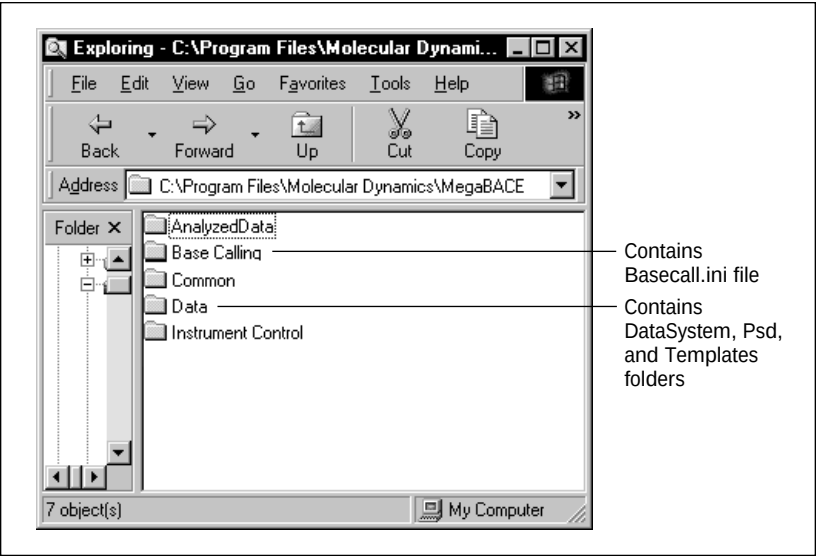


Figure 15-1. Storage locations of the system initialization files and the parameter configuration files.

Table 15-5. The system initialization files

File name	Description	Location
MegaBACE.ini	Defines the options available in the Instrument Control Manager.	...\MegaBACE\Instrument Control\DataSystem folder
Chemistry.ini	Defines the available chemistry parameter sets, including the bases, filters, and laser modes.	...\MegaBACE\Instrument Control\DataSystem folder
Basecall.ini	Defines the default base caller and lists all the base callers available in the Instrument Control Manager.	...\MegaBACE\Base Calling folder

15.4 Parameter configuration file storage

Figure 15-1 shows the system initialization file and parameter configuration file locations. Table 15-6 lists the parameter configuration files that you can use.

Table 15-6. The parameter configuration files

File extension	Description	Storage location
.psd	(Sequencing and STR genotyping only) Used for a given plate, a plate setup data file (.psd) allows you to specify the plate ID, plate setup parameters, and plate- and well-level attributes. (SNP genotyping only) Used for a given plate set, a master.psd file allows you to specify the plate set ID, plate set setup parameters, and the attributes for the first plate injected. Additional, plate-specific .psd files can specify the attributes for each subsequent plate injected during the run.	...\MegaBACE\Instrument Control\Psd folder
.tpl	(Sequencing and STR genotyping only) Used for multiple plates, a plate setup template file (.tpl) allows you to specify the plate setup parameters. The .tpl file uses the same field names as the .psd file for the plate setup parameters. (SNP genotyping only) Used for multiple plate sets, a plate setup template file (.tpl) allows you to specify the plate set setup parameters. The .tpl file uses the same field names as the .psd file for the plate set setup parameters.	...\MegaBACE\Instrument Control\Templates folder
.icp	For use in multiple runs, an instrument control parameter file (.icp) is used as a template to specify a set of parameters that control the instrument.	...\MegaBACE\Instrument Control\Templates folder

Part six

Appendixes

Appendix A Header formats

This appendix provides the standard header format and describes the extended header format used in the MegaBACE raw sample data files (.rsd) and base-called sample data files (.esd). The topics are—

- The standard header format (section A.1)
- Viewing and editing the data in the extended header (section A.2)
- The basic extended header format (section A.3)
- The extended header labels in sequencing files (section A.4)
- The extended header labels in STR genotyping files (section A.5)
- The extended header labels in SNP genotyping files (section A.6)

A.1 The standard header format

Figure A-1 provides the format of the standard header format for the MegaBACE raw sample files (.rsd).

Offset	Size	Type	Variable	
0	24	char	m_tcSampleName[24] ;	// deprecated
24	30	char	m_tcOperator [30] ;	// unused
54	4	float	m_fRunVoltage ;	
58	4	float	m_fRunTime ;	// minutes
62	4	float	m_fInjectionTime ;	
66	4	float	m_fInjectionVoltage ;	
70	4	float	m_fTemperature ;	
74	4	int	m_nLaserIndex ;	// deprecated
78	4	float	m_fPMTVoltage1 ;	
82	4	float	m_fPMTVoltage2 ;	
86	30	char	m_csBeamSplitterA[30] ;	// deprecated
116	30	char	m_csBeamSplitterB[30] ;	// deprecated
146	30	char	m_csFilter1[30] ;	// deprecated
176	30	char	m_csFilter2[30] ;	// deprecated
206	30	char	m_csFilter3[30] ;	// deprecated
236	30	char	m_csFilter4[30] ;	// deprecated
266	4	float	m_fFillTime ;	
270	4	float	m_fFlushTime1 ;	
274	4	float	m_fFlushTime2 ;	
278	4	float	m_fFlushTime3 ;	
282	4	float	m_fPrerunVoltage ;	
286	4	float	m_fRelaxationTime ;	
290	4	float	m_fPrerunTime ;	
294	64	float	m_fSpecSepMx[4][4] ;	
358	40	char	m_csBaseOrder[4][10] ;	// deprecated
398	22	char	m_InstID[22] ;	
420	30	char	m_csParamSetName[30] ;	// deprecated
450	4	DWORD	m_dMagicNumber ;	// 826561349
454	4	DWORD	m_dExtHeaderOffset ;	// new
458	4	DWORD	m_dExtHeaderLen ;	// new
462	4	float	m_fPreinjectionTime ;	// new
466	4	float	m_fPreinjectionVoltage ;	// new
470	4	float	m_fLowPressureTime ;	// deprecated
474	4	float	m_fUserInputTime ;	// new
478	4	float	m_fSleepTime ;	// new
482	4	float	m_fSleepTemperature ;	// new
486	10	char	m_csBuildNumber[10] ;	//format: 19990927A, includes \0
496	4	int	m_nNumColumns ;	// maximum number of letter labels on the plate (depends on instrument type)
500	4	int	m_nNumRows ;	// maximum number of number labels on the plate (depends on instrument type)
504	4	float	m_fIntervalVoltage	
508	2	WORD	m_nIntervalTime	
510	4	float	m_fPMTVoltage3	
514	4	float	m_fPMTVoltage4	
518	2	WORD	m_wNumberOfCapillaries	
520	64	BYTE	m_byReserved[48] ;	
568	6	char	m_csRunID[6] ;	
574	18	char	m_csPlateID[18] ;	// deprecated
592	4	DWORD	m_dwPrevRunID ;	// unused
596	4	BOOL	m_bDWORDFile ;	// 1=raw,0=processed
600	4	BOOL	m_bSaturated ;	
604	4	float	m_fVersion ;	// 2.000
608	4	time_t	m_tRunStartTime ;	
612	4	time_t	m_tRunStopTime ;	
616	4	DWORD	m_dwNumberOfLines ;	
620	4	int	m_nScanRate ;	// 0=1.75, 1=3.0
624	1	BYTE	m_byFileType ;	// reserved
625	1	BYTE	m_byNumberOfCapillaries ;	// deprecated
626	1	BYTE	m_byCapUsage ;	
627	1	BYTE	m_byDummy ;	
628	4	UINT	m_uiSecurityCode ;	

Figure A-1. Standard header format.

A.2 Viewing and editing the data in the extended header

You use the MegaBACE Header Editor software to open the file and display the contents of the extended header. If necessary, you can edit information, such as the chemistry names.

To start MegaBACE Header Editor, in the **MegaBACE** shortcut folder on the Windows desktop, double-click the **Header Editor** icon. The Header Editor window appears (figure A-2).

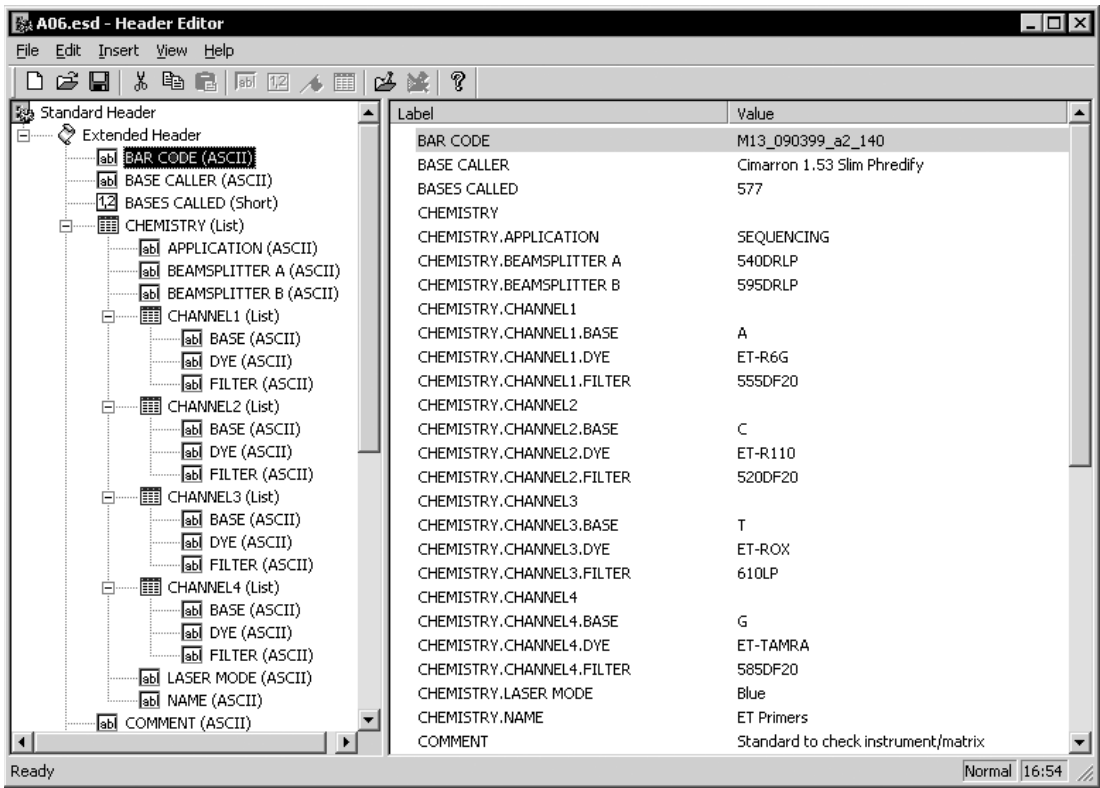


Figure A-2. The Header Editor window displaying the information for a selected well from a sequencing run.

The MegaBACE Header Editor software includes Help that describes how to use the software. To access Help, in the MegaBACE Header Editor software, choose **Help Topics** on the **Help** menu.

A.3 The basic extended header format

The extended header is a variable length, variable content data block that is encoded using the Amersham Biosciences self-describing data protocol (appendix B).

All entries in the extended header are label-value pairs. The label is always an ASCII string. The values may be in any supported format.

The value associated with a label may be an array or a list. The entries in a list are also packed as label-value pairs. Therefore, the length of any list is always double the number of entries in the list.

The entries in the extended header are alphabetized by label. The entries within a list are also alphabetized.

The application programming interface (API) for retrieving values from the extended header accesses items in a list using a dot-delimited hierarchical name.

For example, CHEMISTRY.CHANNEL1.BASE would return the value associated with the label BASE from the CHANNEL1 list inside the CHEMISTRY value. (Labels are treated as case-independent.)

Any user-defined data imported from a plate setup data file (.psd) will be included in the extended header.

A.4 The extended header labels in sequencing files

For sequencing runs, you can view the extended header labels in the raw sample data files (.rsd) and the base-called data files.

A.4.1 The extended header in the sequencing .rsd files

The following labels will always appear in the raw sample data files (.rsd) for sequencing runs.

Note: The indentations within CHEMISTRY denote nested lists.

BAR CODE
BASE CALLER
CHEMISTRY
 APPLICATION
 BEAMSPLITTER A
 BEAMSPLITTER B
 BEAMSPLITTER C
 CHANNEL1
 BASE
 DYE
 FILTER
 CHANNEL2
 BASE
 DYE
 FILTER
 CHANNEL3
 BASE
 DYE
 FILTER
 CHANNEL4
 BASE
 DYE
 FILTER
 LASER MODE
 NAME
COMMENT
MACHINE ID
PLATE ID
SAMPLE NAME
WELL ID

A.4.2 The bar code

BAR CODE is the string entered into the Plate ID field on the Plate Setup window. The bar code can be entered either by keyboard entry or with a bar code scanner. If there is no .psd file associated with the name, or if there is no PLATE ID entry in the .psd file, the plate ID will match the bar code.

A.4.3 The extended header in base-called sample data files

All entries in a raw sample data file (.rsd) are preserved in the corresponding base-called sample data file (.esd). The following additional entries appear in base-called sample data files:

BASES CALLED
END GOOD SEQUENCE
END SEQUENCE
END SIGNAL
END TRACE PROCESSING
PEAK POSITIONS
QUALITY INDEX
QUALITY SCORES
SEQUENCE
START GOOD SEQUENCE
START SEQUENCE
START SIGNAL
START TRACE PROCESSING

A.5 The extended header labels in STR genotyping files

The following labels will always appear in the raw sample data files (.rsd) for STR genotyping runs.

Note: The indentations within CHEMISTRY denote nested lists.

BAR CODE
CHEMISTRY
 APPLICATION
 BEAMSPLITTER A
 BEAMSPLITTER B
 BEAMSPLITTER C
 CHANNEL1
 DYE
 FILTER
 LASER MODE
 NAME
COMMENT
MACHINE ID
PLATE ID
SAMPLE NAME
WELL ID

Note that BAR CODE is the string entered into the Plate ID field on the Plate Setup screen, either by keyboard entry or with a bar code scanner. If there is no .psd file associated with the name, or if there is no PLATE ID entry in the .psd file, the plate ID will match the bar code.

A.6 The extended header labels in SNP genotyping files

The following labels will always appear in the raw sample data files (.rsd) for sequencing runs.

A.6.1 Additional extended header labels for SNP genotyping

The labels that appear in the raw sample data files (.rsd) from a SNP genotyping run can include—

- **Plate set attributes**—The attributes that apply to the entire run appear at the top level of the extended file header. Table A-1 lists the additional top-level labels for SNP genotyping.

- **Injection-specific (plate) attributes**—The attributes for an injection appear under the subsection for the given injection. Injection-specific subsections are labeled as INJECTION 1, INJECTION 2, and so forth. Table A-2 lists the labels for the injection-specific attributes.

Table A-1. Plate set fields in the extended header

Field	Description
CHEMISTRY APPLICATION BEAMSPLITTER A BEAMSPLITTER B BEAMSPLITTER C CHANNEL BASE DYE FILTER LASER MODE NAME	Always present. The chemistry parameters are a plate set definition setting (section 12.2). Note: The indentations within CHEMISTRY denote nested lists.
COMMENT	(Optional) Plate set definition setting (section 12.2).
MACHINE ID	Always present, the software adds the ID of the MegaBACE instrument.
WELL ID	Always present, the well ID associated with the corresponding .rsd file.
INTERVAL VOLTAGE	Always present, a plate set definition setting (section 12.2).
INTERVAL TIME	Always present, a plate set definition setting (section 12.2).
NUMBER OF INJECTIONS	Always present, the actual number of sample injections performed during the run.
PLATE SET ID	Always present. If a .psd file is used for the plate set and the plate set ID is missing, the software uses the plate ID, if available. Otherwise, the software uses the bar code or the .psd file name.
PLATE SET	(Optional) Values may be present if a .psd file is used (section 13.2.2).
All user-defined plate set labels	(Optional) Values may be present if a .psd file is used.

Table A-2. Injection-specific fields in the extended header

Field	Description
BAR CODE	Injection-specific attribute for SNP genotyping, but a top-level, run-wide attribute in other applications, such as sequencing and STR genotyping. The value may be missing for injections other than the first injection if a .psd file is not used.
PLATE ID	(Optional) Injection-specific attribute for SNP genotyping, but a top-level, run-wide attribute in other applications, such as sequencing and STR genotyping. The value may be missing.
MARKER NAME	Always present for each injection of a SNP genotyping run. If only one SNP marker name is provided, the software copies that name for each injection in the run.
SAMPLE NAME	(Optional) Injection-specific attribute for SNP genotyping, but a top-level, run-wide attribute in other applications, such as sequencing and STR genotyping. For SNP genotyping, a .psd file must be used to specify sample names.
All user-defined labels	(Optional) Values may be present if a .psd file is used.

Appendix B Protocol for self-describing data

This appendix describes a protocol for packing self-describing data into a file. The topics are—

- Introduction (section B.1)
- The 1-byte header (section B.2)
- The length value (section B.3)
- The data buffer (section B.4)
- Notes about the protocol (section B.5)

B.1 Introduction

The extended header consists of a list of keyword-value pairs. However, the protocol is completely general and can be used to pack data of almost any form.

In the protocol, a data block consists of a series of elements, where each element contains:

- A 1-byte header that contains two fields.
- A length value that consists of 1 to 3 bytes (except for Terminator, which is always zero) as determined by the number of length bytes field in the header.
- A data buffer that contains the number of bytes specified in the length value.

B.2 The 1-byte header

The 1-byte header contains two fields:

- A 6-bit format field that defines the type of data contained in the element. The format field consists of the most significant 6 bits of the header byte.
- A 2-bit field that states the number of length bytes used to specify the size of the data storage. The number of length bytes field consists of the least significant 2 bits of the header byte. This field is zero for Terminator. For all other formats it must be nonzero.

Figure B-1 shows the protocol.

Format (octal)	Symbol	Description
00	Ln	List (where n is a number)
01	An	ASCII
02	Cn	Unicode
10	B	Boolean (1 byte: 1=true, 0=false)
31	U1	1-byte unsigned integer (unsigned char)
[32	U2	2-byte unsigned integer (unsigned short)]
[34	U4	4-byte unsigned integer (unsigned long)]
40	F8	8-byte floating (double)
44	F4	4-byte floating (float)
[51	S1	1-byte signed integer (char)]
52	S2	2-byte signed integer (short)
54	S4	4-byte signed integer (long)
77	T	Terminator

Figure B-1. The protocol for self-describing data.

Notes to figure B-1: The symbol that represents a format is used only in readable dumps of the data block. The symbol is not embedded in the data structure itself.

The lines delimited with square brackets indicate format codes that are not allowed in the extended header because those data types are not supported by the (COM) interface that enables the header.

B.3 The length value

The length value consists of 1 to 3 bytes (except for Terminator) as determined by the number of length bytes field in the header.

If the format is List, this field contains the number of elements in the list. Otherwise, it contains the number of bytes of data storage used in the element.

Note: The number of individual data values in the element (that is, the array size) is the length divided by the size of each value.

B.4 The data buffer

For all formats except List, the data buffer contains the number of bytes specified in length. For List, the number of elements is specified in length.

B.5 Notes about the protocol

The protocol packs the elements as a byte stream. In all cases, multibyte data are packed with the least significant byte first. This means that the packing/unpacking software will differ depending on whether the platform implements a big endian or little endian architecture (x86 is little endian).

The size of the data block may not be known a priori. If the size is indeterminate, the Terminator format code can be used to tell the unpacking algorithm where to stop. If the size of the data block can be unambiguously determined, the use of the Terminator is optional.

Character strings (A and C formats) are null terminated. (See figure B-1 for notations.)

The List format is used for creating complex data structures, which are treated as single elements. If the value requires multiple parameters of different data types (for example, character, string, or int, with int to represent the base, filter center, and filter width for each channel in the Chemistry specifier), the parameters can be specified inside a List, thereby making them a single value associated with the corresponding Keyword.

Human-readable representation (dump format) is as follows:

- By convention, the format and data for each element are listed between angle brackets. If the element is represented on a single line, the brackets appear on the same line. Otherwise, the brackets are vertically aligned. Nesting is indicated by indentation.
- For List elements, the number of elements in the List is appended to the L format identified. For example, L3 denotes a list containing three elements (figure B-2).

Note: An embedded List counts as a single element for its enclosing List. Lists may be nested to any depth.

```
< L3
    < A6   Gizmo >
    < S2   -414 >
    < L2
        < ... >
        < ... >
    >
>
```

Figure B-2. Example of the List format. (See figure B-1 for notations.)



For character elements (A and C formats), the number of characters is appended to the format identifier. For example, A6 denotes a five-character ASCII string with its terminating null. Strings are not surrounded by quotes unless the quote characters are supposed to be encoded as part of the string.

All other data types are treated as arrays if the element contains more than one item, where the array size is shown in square brackets. For example, U2[4] indicates an array of four 2-byte unsigned integers (figure B-3).

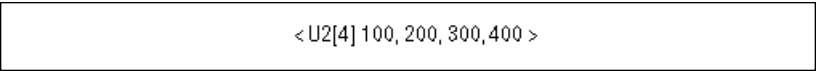


Figure B-3. Example of an array. (See figure B-1 for notations.)



Appendix C Launching a cleanup script after a run

Important

To write scripts that are launched using the CLEANUP COMMAND field, you should be familiar with Windows Shell programming.

The CLEANUP COMMAND field in the MegaBACE.ini file allows you to specify a program or batch file to be executed at the end of the run for the selected application (sequencing, STR genotyping, or SNP genotyping). (For the sequencing application, the program or batch file can also be executed after automated base calling.) For instance, you can use the CLEANUP COMMAND feature to transfer plate records and raw data files to another workstation for analysis or archiving. This appendix describes how to use the CLEANUP COMMAND field in the MegaBACE.ini file. The topics are—

- Using the CLEANUP COMMAND (section C.1)
- Batch file examples (section C.2)

C.1 Using the CLEANUP COMMAND

To use the CLEANUP COMMAND feature:

1. Create the program executable file (.exe) or batch file (.bat) that you want the Instrument Control Manager to launch. See section C.2 for examples of batch files (.bat).
 - Make use of the program arguments that the Instrument Control Manager supplies at run time. See tables C-1 and C-2 for a list of arguments.
 - Make a note of the location where you store the .exe or .bat file.
2. Open the MegaBACE.ini file, and locate the section for the application you are running (Sequencing, Genotyping, or SNP). In the value column for the CLEANUP COMMAND field, type the file path of the .exe or .bat file. For example:

Field in MegaBACE.ini file	Value
CLEANUP COMMAND	..\test.bat

Make sure you use a tab (or tabs) between the CLEANUP COMMAND field and the text you type.



The file path that you type in the Value column must be the full path of the .exe or .bat file or a path relative to the MegaBACE\Instrument Control folder. For example, if the file is stored in:

- **"...\MegaBACE\Instrument Control" folder**—Type only the file name, for example, test.bat.
- **Any folder other than the Instrument Control folder**—Type either the:

Relative path, for example, ..\test.bat for a file stored in the MegaBACE folder

Full path, for example, D:\test.bat

Note: You do not necessarily need to use quotes around the path even if it contains spaces.

3. Restart the Instrument Control Manager, and perform the run. See the *MegaBACE 4000 Instrument User's Guide* for details.

At the end of the run (after base calling, if any, for sequencing), the Instrument Control Manager runs the command with the program arguments (tables C-1 and C-2), which the Instrument Control Manager provides dynamically. The command runs in the background (without displaying a window).

Table C-1. Instrument Control Manager arguments for the sequencing and STR genotyping applications

Argument number	Description*
1	The plate ID.
2	The full path of the run folder (in quotes).
3	The application name, Sequencing or Genotyping. Note: The application name appears in the title bar of the Instrument Control Manager windows.
4	The bar code of the plate. This is the text that was entered in the Plate ID field during plate setup.

* Items allowed to contain spaces are passed in quotes.



Table C-2. Instrument Control Manager arguments for the SNP genotyping application

Argument number	Description*
1	The plate set ID.
2	The full path of the run folder (in quotes).
3	The application name, SNP Genotyping. Note: The application name appears in the title bar of the Instrument Control Manager windows.
4	The bar code used for the plate set definition. This is the text that was entered in the Plate Set ID field during plate set setup.

* Items allowed to contain spaces are passed in quotes.

C.2 Batch file examples

This section describes some simple examples of batch files (.bat) that could be specified for the CLEANUP COMMAND field.

C.2.1 Displaying the full set of arguments after the run

In the test.bat file example (figure C-1a), the code specifies that the full set of arguments, supplied by the Instrument Control Manager, will appear on screen in a Notepad file entitled, a.txt. Figure C-1b shows the resulting file that appears after the run.

In order for the Instrument Control Manager to launch the test.bat file (figure C-1a) after the run, the .bat file name and file path must be specified in the MegaBACE.ini file (section C.1), and the .bat file must reside at the specified location.

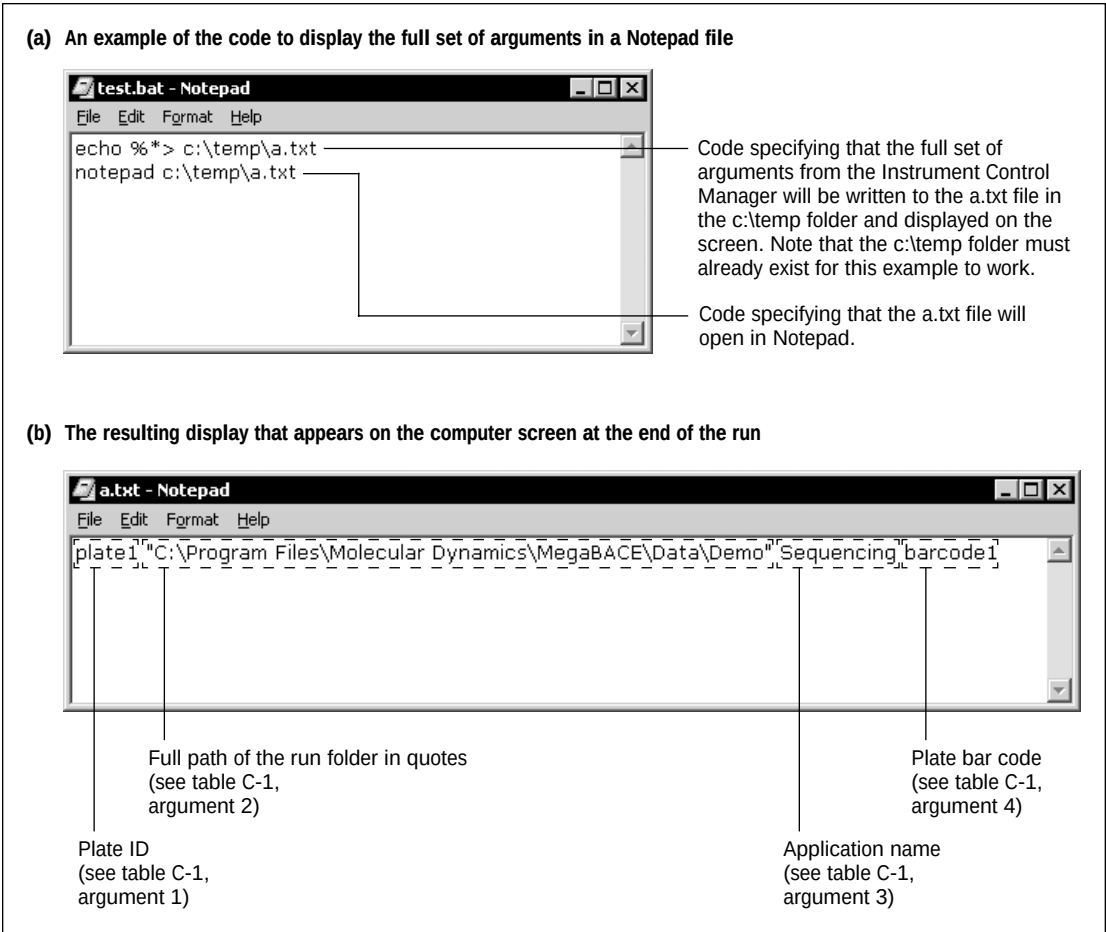


Figure C-1. An example of a .bat file displayed in Notepad (a), and the on-screen display that appears as a result of the code (b).

C.2.2 Coping the raw sample data files to another folder

Figure C-2 shows an example of the code to copy the raw sample data files (.rsd) from the run folder to another folder after the run. Figure C-2 shows the .bat file name and file path specified in the MegaBACE.ini file.

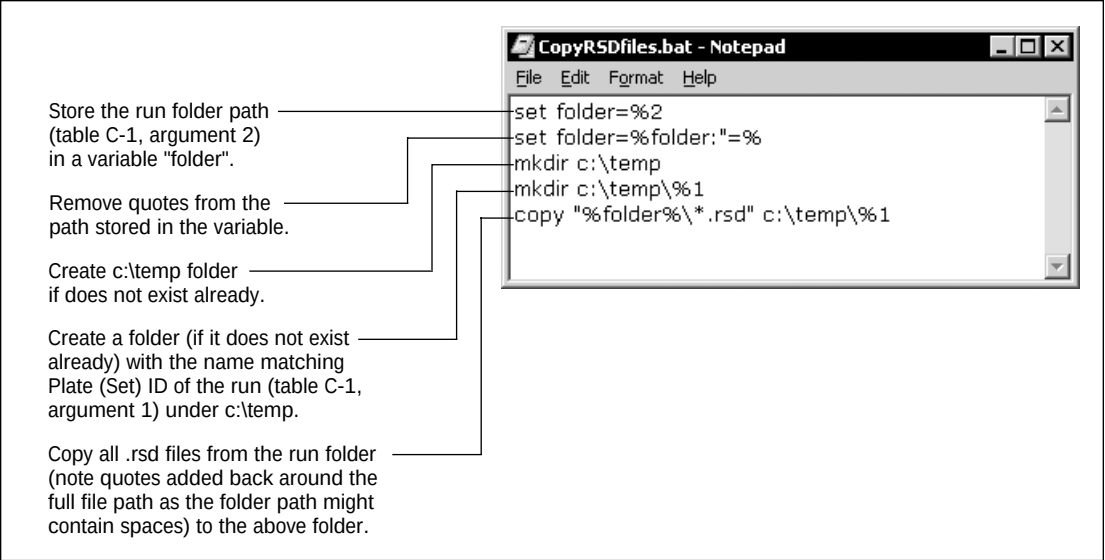


Figure C-2. An example of a .bat file displayed in Notepad and the functional explanation for each line of code.

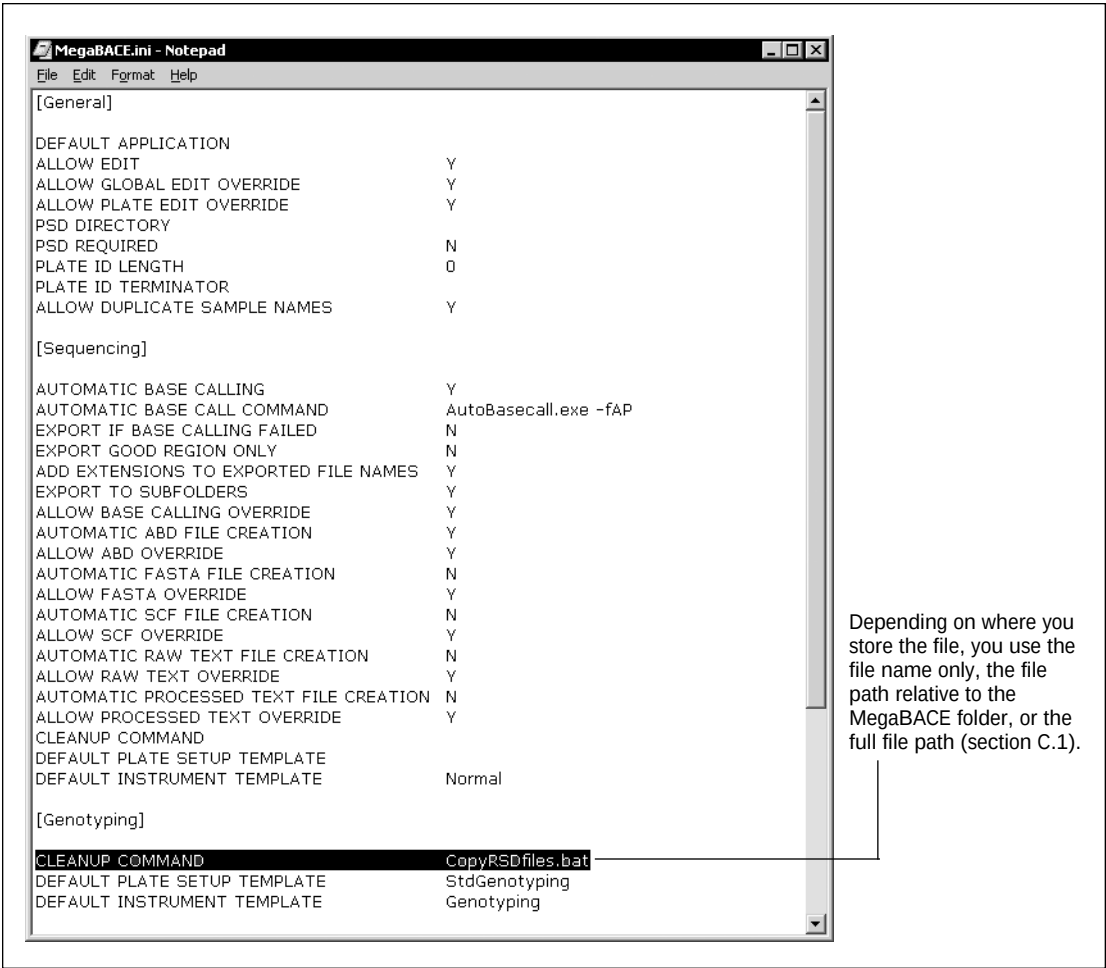


Figure C-3. The MegaBACE.ini file showing the CLEANUP COMMAND field for STR genotyping. Note that only the file name is required if the .bat file is stored in the MegaBACE\Instrument Control folder.

Glossary

.abd files—(sequencing only) the base-called data files that can be viewed using ABD software and used by the Phred application. Each .abd file contains the data for a given well on a plate (for example, A01), including the raw and analyzed electropherogram data.

.esd files—(sequencing only) the base-called data files. Each .esd file contains the called sequence, read length, sequencing starting and ending points, quality values, current profile, and the analyzed electropherogram for a given well on a plate (for example, A01). The sequenced data can be generated automatically by the Instrument Control Manager after a sequencing run or by the Sequence Analyzer software. Note that the Instrument Control Manager stores the .esd files for each run in an analyzed run folder.

.icp files—the instrument control template files, which can be used as templates to specify the instrument run conditions and the matrix fill and flush cycles for multiple runs.

.psd files—the plate setup data files that you can use to automatically import attributes for a plate and for individual wells on the plate, such as sample names. The Instrument Control Manager includes the information from the .psd file in the extended header of each raw sample data file (.rsd).

(SNP genotyping only) The plate setup data files allow you to specify attributes for a plate set or a plate.

.rsd files—the raw sample data files, each of which contains the raw data for a given well on a plate (for example, A01). The storage location is in a corresponding raw run folder in the ...\\MegaBACE\\Data folder (default) or the location you specify.

.scf files—(sequencing only) the base-called data files in SCF format. Each .scf file contains the raw data and electropherogram. If analyzed, each .scf file includes the raw and base-called data and the called sequence.

.seq files—(sequencing only) the base-called data files in FASTA format. Each .seq file contains the .esd file name, .esd file location, and the called sequence for a given well on a plate. You can view these files in another application, such as Notepad.

Glossary

.tpl files—the plate setup template files. For sequencing and STR genotyping, the .tpl file can specify the plate setup parameters, including the electrophoresis conditions, the chemistry parameters, the file names, and optional parameters, such as the base caller.

For SNP genotyping, the .tpl file can specify the plate set setup parameters, which apply to the entire set of plates injected during a single run.

analyzed run folder—(sequencing only) the folder that contains the base-called sample files for a given run. The Instrument Control Manager uses the plate ID and base caller ID to name the analyzed run folder. The Instrument Control Manager creates an analyzed run folder for the data from each plate for which the Instrument Control Manager performs automatic base calling. The software stores the resulting analyzed run folder in the ...\\MegaBACE\\AnalyzedData folder (default) or the location you specify.

base callers—(sequencing only) the software that identifies the candidate peaks and calls the bases in a sequence.

base colors—the display color of the base letter and trace associated with a given base: A = green, C = blue, T = red, G = black. Although the conventional color representation for G is yellow, the trace and letter are displayed in black for ease of viewing.

channels—see spectral channel.

chemistry parameters—the names of the dyes, the base order or dye-to-channel mapping, the names of the filters and beamsplitters, and the laser mode used for a given application (sequencing, STR genotyping, or SNP genotyping).

dye set—the fluorescent dyes used to label your experiment. The MegaBACE system is capable of detecting and separating the emissions of multiple dyes of different colors in a single capillary.

electropherograms—the digitized graphs the system produces from the fluorescent intensity detected from the capillaries during the run. The system produces an electropherogram for each capillary. Each electropherogram consists of four colored traces that represent the signals detected through the four spectral channels.

electrophoresis parameters—the settings for the sample injection voltage and time and the run voltage and time. The electrophoresis parameters are displayed in the Electrophoresis tab in the Plate Setup or Plate Set Setup window of the Instrument Control Manager software.

instrument parameters—a combination of settings defining the instrument run conditions and matrix fill and flush cycles. The instrument parameters are displayed in the Instrument Control window of the Instrument Control Manager software.

marker—see SNP marker of STR marker.

multiple injection—a proprietary method of increasing SNP detection throughput by a more efficient use of capillary volume. Multiple plates are injected into the same set of capillaries, at time-spaced intervals, during the waiting period before the first injected samples appear at the detector.

Phred—(sequencing only) sequence read editor program used to verify the accuracy of sequenced DNA. The program was developed by Phil Green at the University of Washington.

plate definition—includes the plate ID and plate or plate set setup parameters.

plate ID—the name you give to the plate when you create a plate definition in the Instrument Control Manager. For SNP genotyping, the plate ID is optional.

plate set—(SNP genotyping only) a group of plates that share plate setup parameters and are serially injected into the same set of capillaries during a single run on the instrument.

plate set definition—(SNP genotyping only) the combination of the plate set ID and plate set setup parameters for a group of plates that you inject during a multi-injection run. The plate set definition can also include the list of bar codes for plates in the plate set and the injection order.

plate set ID—(SNP genotyping only) the name you give to the set of plates you inject in a multi-injection run. The software uses the plate set ID to name the raw run folder that stores the data from the run. If you do not provide a plate set ID, the software uses the ID of the plate that is injected first as the ID for the entire plate set. If the plate ID is also missing, the software uses the master .psd file name (bar code).

plate set setup parameters—(SNP genotyping only) the plate setup parameters for a multiple-injection run, where all the setup parameters are common to all injections (plates). The plate set setup parameters can include electrophoresis conditions, the chemistry parameters, the file names, and optional parameters. In addition, the parameters can include the list of all the plates (plate IDs or plate bar codes) to be injected during the run.

plate setup parameters—a combination of electrophoresis parameters, chemistry parameters, sample names, and optional parameters that define a plate.

Glossary

raw data—the data collected by the instrument that have not been analyzed. The Instrument Control Manager software creates a raw run folder for the raw sample data files (.rsd) for each plate you run.

raw run folder—the folder that contains the raw sample data files (.rsd) for a given run. The Instrument Control Manager stores the raw run folder containing the associated .rsd files in the ...\\MegaBACE\\Data folder (default) or the location you specify.

run—the process of injecting and scanning a plate of samples on the MegaBACE instrument and detecting the resulting signal from each capillary. Each run has a unique date and user ID.

run folder—see analyzed run folder or raw run folder.

run ID—a unique designation the Instrument Control Manager software assigns to each run of a sample plate on the instrument.

sample files—see .rsd files. For sequencing, see also .abd files, .esd files, .seq files, and .scf files.

sample name—(sequencing and STR genotyping only) the designation you provide for a given well in the Sample Names tab or in a plate setup template or .psd file.

(SNP genotyping only) the sample name is the designation you provide for a given well for a given injection (plate) in the multi-injection run. A plate setup data (.psd) file is required for each injection for which you want to assign sample names.

scan number—a number representing a sampling of the data during the run and describing the location of a data point. Instrument Control Manager samples the data continuously during a run at a rate of 1.75 Hz (105 times per minute).

SNP (single nucleotide polymorphism)—(SNP genotyping only) a base position in the genome that is variable in a population. To detect SNPs on the MegaBACE instrument, you use the SNUPe genotyping kit to assay your samples, you perform a multi-injection run, and then you use SNP Profiler to perform the SNP genotyping analysis.

SNP marker—(SNP genotyping only) any genetic locus containing a single-base variation or polymorphism.

SNUPe genotyping kit—(SNP genotyping only) a set of components used to assay SNP samples based on single-nucleotide primer extension (SNUPe) technology.



spectral channel—the combination of laser, beamsplitter, optical filter, and PMT the system uses to detect the emission signals of a given dye. The MegaBACE instrument has a total of four spectral channels. The output of each spectral channel is represented by a different colored trace in the electropherogram.

STR marker—short-tandem-repeat (STR) marker, also known as a microsatellite marker, refers to any genetic locus containing a repeating variation on a chromosome. An STR marker consists of two primers that are used in the PCR reaction.

traces—the four curves in the electropherogram representing the signals detected through the four spectral channels of the instrument. Each trace is displayed in a specific color to represent the associated base or dye.

Index

.abd files. *see* ABD files
.esd files. *see* base-called sample data files
.icp files. *see* instrument control templates (.icp)
.psd files. *see* plate setup data files (.psd)
.rsd files. *see* raw sample data files (.rsd)
.scf files. *see* SCF files
.seq files. *see* FASTA files
.tpl files. *see* plate setup templates (.tpl)
or plate set setup templates (.tpl)

A

ABD files
 MegaBACE.ini fields 5-4
 storage 15-4
administrator tasks
 sequencing 1-3
 SNP genotyping 1-13
 STR genotyping 1-8
Allow Edit field 2-4
Allow Global Edit Override field 2-4
Allow Plate Edit Override field 2-4
annotation information for base calling 5-12
applications
 see also sequencing, STR genotyping, or SNP genotyping
 about changing 1-2
 default, specifying 2-4, 2-5
assistance xiv

B

bad capillaries, specifying in MegaBACE.ini file 2-10 to 2-11
bar codes
 automating bar code entry 2-9
 Plate ID Length and Plate ID Terminator fields 2-4
 sequencing
 extended file header A-6
 importing a .psd file 7-2
 SNP genotyping
 extended file header A-9
 importing a .psd file 13-4
 PLATE SET field 11-9

Index

STR genotyping
 extended file header A-7
 importing a .psd file 10-2
base call command, automatic 5-3
base callers
 methods to specify 5-12 to 5-13
 specifying in .psd file 7-7 to 7-8
base calling, automatic
 annotation information 5-12
 configuring MegaBACE.ini 2-4, 5-7 to 5-10
 file storage 15-3 to 15-5
 overview 1-4
 setting up 5-7 to 5-13
 turning on or off 5-7
Basecall.ini file
 fields 5-12
 format overview 1-20
 sections and fields 5-11 to 5-12
 storage 1-20, 15-6
 using 5-8 to 5-12
base-called sample data files
 creation 5-9 to 5-10
 storage 15-4

C

capillaries, empty
 specifying in .psd file
 sequencing 7-10
 SNP genotyping 13-9
 STR genotyping 10-8
 specifying in plate set setup parameters
 (SNP genotyping) 12-7 to 12-8
 specifying in plate setup parameters
 sequencing 6-7 to 6-8
 STR genotyping 9-6 to 9-7
changing applications, about 1-2
Chemistry Parameters tab
 sequencing 6-4
 SNP genotyping 12-4
 STR genotyping 9-4

Chemistry.ini file
 adding new chemistry parameters 3-3 to 3-5
 configuring 3-1
 errors 3-5
 fields 3-2 to 3-3
 format overview 1-20
 storage 1-20, 15-6
Cleanup command
 sequencing 5-5
 SNP genotyping 11-3
 STR genotyping 8-3
 using the command C-1 to C-6
COM interface B-2
Comments tab
 sequencing 6-5
 STR genotyping 9-4

D

data storage
 sequencing
 base called data 15-3 to 15-5
 raw sample data files (.rsd) 15-1 to 15-2
 SNP genotyping 15-2 to 15-3
 STR genotyping 15-1 to 15-2
defaults
 application, specifying 2-4, 2-5
 sequencing
 base caller, setting 5-8 to 5-12
 instrument control template 5-5, 5-6
 plate setup template 5-5 to 5-6, 6-8 to 6-9
 SNP genotyping
 instrument control template 11-7
 plate set setup template 11-6, 12-9
 STR genotyping
 instrument control template 8-3, 8-4 to 8-5
 plate setup template 8-3 to 8-4, 9-8
documentation, user xiii
duplicate sample names, allowing 2-5

E

edit mode

- allow edit 2-4
- Allow Global Edit Override field 2-4
- Allow Plate Edit Override field 2-4
- enabling or disabling 2-6 to 2-7

Electrophoresis Parameters tab

- sequencing 6-3, 6-4
- SNP genotyping 12-3, 12-4
- STR genotyping 9-3, 9-4

empty capillaries

- specifying in .psd file
 - sequencing 7-10
 - SNP genotyping 13-9
 - STR genotyping 10-8
- specifying in plate set setup parameters (SNP genotyping) 12-7 to 12-8
- specifying in plate setup parameters
 - sequencing 6-7 to 6-8
 - STR genotyping 9-6 to 9-7

Enter Password window 14-2

export files, storage locations 15-3, 15-4

export options, sequencing 5-9 to 5-10

extended file header

- basic format A-4
- sequencing files A-4 to A-6
- SNP genotyping files A-7 to A-9
- STR genotyping files A-7

F

FASTA files

- MegaBACE.ini files 5-4
- options 5-10
- storage 15-4

features, Instrument Control Manager 1-2

file export, sequencing 5-9 to 5-11

file names

- .psd files 7-1, 10-1, 13-4
- sequencing 7-9, 15-1 to 15-2, 15-3 to 15-5

- SNP genotyping 13-2, 13-6, 15-2 to 15-3
- STR genotyping 10-7, 15-1 to 15-2

File Names tab (SNP genotyping) 12-4

files, initialization

- see *also* Chemistry.ini file, Basecall.ini file, or MegaBACE.ini file
- format overview 1-20
- storage 1-20, 15-6

files, parameter configuration

- see *also* instrument control templates, plate setup or plate set setup templates, or plate setup data files
- format overview 1-20
- storage 1-20, 15-7

files, sample data

- protocol for self-describing data B-1 to B-4
- sequencing
 - .rsd file storage 15-1 to 15-2
 - base-called file storage 15-3 to 15-5
 - extended header A-4, A-4 to A-6
 - file creation options 5-4 to 5-5
 - format options 5-9
- SNP genotyping
 - extended header A-4, A-7 to A-9
 - naming and storage 15-2 to 15-3
- standard header A-1
- STR genotyping
 - extended header A-4, A-7
 - naming and storage 15-1 to 15-2

files, suppressing creation

- bad capillaries 2-10 to 2-11
- empty capillaries
 - sequencing 6-7 to 6-8, 7-10
 - SNP genotyping 13-9
 - STR genotyping 9-6 to 9-7, 10-8

folders

- storage 15-1 to 15-7
- subfolders for sequencing, specifying 5-10

Index

H

Header Editor 1-1
Header Editor software A-3
header, file
 extended A-4 to A-9
 standard A-1
Host Scan Controller 1-1

I

importing .psd files
 automating bar code entry 2-9 to 2-10
 MegaBACE.ini fields 2-4
 requiring .psd files 2-8
 storage location, specifying 2-8
initialization files. see Basecall.ini file, Chemistry.ini file, or MegaBACE.ini file
Instrument Control Manager software 1-1, 1-2
instrument control parameters
 creating templates 4-1 to 4-6
 fields 4-3
instrument control templates (.icp)
 creating from a text file 4-1 to 4-4
 creating in the Instrument Control window 4-1
 setting default
 sequencing 5-5, 5-6
 SNP genotyping 11-3, 11-7
 STR genotyping 8-3, 8-4
 storage 15-7

L

List of Protocols area 14-1

M

MegaBACE.ini file
 about sections in file 2-1
 bad capillaries, specifying 2-10 to 2-11
 example 2-2
 format overview 1-20

general section

 application, specifying default 2-5 to 2-6
 edit mode, enabling or disabling 2-6 to 2-7
 field names and data types 2-4 to 2-5
 modifying 2-3

sequencing

 automatic base calling and file export 5-7 to 5-10
 fields 5-1 to 5-5
 fields, illustrated 5-2

SNP genotyping

 fields 11-1 to 11-5
 fields, illustrated 11-2
 injection specifications 11-7 to 11-9
 plate information, prompting operator 11-9 to 11-13

storage 1-20, 15-6

STR genotyping

 fields 8-1 to 8-3
 fields, illustrated 8-2

multi-injection marker 1-15

O

Optional Parameters tab

 sequencing 6-4, 6-5
 SNP genotyping 12-5
 STR genotyping 9-4

Override Commands window 14-2

P

plate attributes (SNP genotyping)

 .psd files 13-3 to 13-5, 13-6, 13-10 to 13-11
 overview 1-18
 prompting the operator 11-9 to 11-13

plate definitions

sequencing

 creating, methods for 6-1 to 6-2
 workflow alternatives overview 1-5 to 1-6

STR genotyping

 creating, methods for 9-1 to 9-2
 workflow alternatives overview 1-10 to 1-11

- plate IDs
 - bar codes, MegaBACE.ini file settings for 2-9
 - MegaBACE.ini file settings 2-4
 - sequencing 7-1
 - SNP genotyping
 - prompting operator 11-4, 11-12
 - specifying in .psd files 13-6, 13-10
 - STR genotyping 10-1
- plate set definitions (SNP genotyping)
 - creating master .psd files 13-1 to 13-2, 13-4 to 13-9
 - creating, methods for 12-1 to 12-2
 - defining plate list in master .psd file 13-8
 - workflow alternatives overview 1-15 to 1-16
- plate set IDs (SNP genotyping) 13-7
- plate set setup parameters (SNP genotyping)
 - fields in .psd files 13-5 to 13-6
 - overview 1-17
 - plate set setup templates 12-1 to 12-9
- plate set setup templates (SNP genotyping)
 - creating in Plate Set Setup window 12-5 to 12-6
 - creating in text file 12-6 to 12-7
 - creating plate set definitions, about 12-2
 - default, setting 11-3, 11-6, 12-9
 - described 12-4 to 12-5
 - illustrated 12-7
 - storage 15-7
- Plate Set Setup window (SNP genotyping)
 - creating plate set definitions 12-2
 - illustrated 12-3
- plate setup data files (.psd)
 - automating import 2-9 to 2-10
 - fields in MegaBACE.ini file 2-4
 - format overview 1-20
 - requiring 2-8
 - sequencing
 - base callers, specifying 7-7 to 7-8
 - creating 7-4 to 7-6
 - described 7-1
 - field descriptions 6-4, 7-6
 - file name requirements 7-1 to 7-3
 - importing, about 6-2
 - workflow overview 1-5
 - SNP genotyping
 - creating 13-4 to 13-5
 - described 13-1
 - importing, about 12-2
 - master file for plate set 13-1, 13-2, 13-4 to 13-9
 - plate attributes 13-3, 13-4, 13-6, 13-10
 - storage 1-20, 15-7
 - storage location, specifying 2-8
 - STR genotyping
 - creating 10-4 to 10-6
 - described 10-1
 - field descriptions 9-4
 - fields 10-6
 - file name requirements 10-1 to 10-3
 - importing, about 9-2
- plate setup parameters
 - sequencing
 - default template, setting 5-5 to 5-6
 - defined 6-4
 - fields in .psd files 7-6
 - STR genotyping
 - default template, setting 8-3 to 8-4
 - defined 9-4
 - described 9-3 to 9-4
 - fields in .psd files 10-6
- plate setup templates (.tpl)
 - sequencing
 - creating in text file 6-6
 - creating plate definitions, about 6-2
 - default, setting 6-8 to 6-9
 - described 6-3 to 6-5
 - storage 15-7
 - STR genotyping
 - creating in Plate Setup window 9-5
 - creating in text file 9-6
 - creating plate definitions 9-2
 - default, setting 9-8
- Plate Setup window
 - sequencing 6-2, 6-3
 - STR genotyping 9-2, 9-3

Index

protocol for self-describing data B-1 to B-4
protocols, instrument
 bypassing a protocol 14-2
 List of Protocols area 14-1
 overriding a protocol 14-3
publications xiii

R

raw sample data files (.rsd)
 sequencing 15-1 to 15-2
 SNP genotyping 15-2 to 15-3
 STR genotyping 15-1 to 15-2

S

Sample Injection window, illustrated 11-8
sample names
 sequencing
 .psd file 7-9
 plate setup template 6-4
 Plate Setup window 6-5
 SNP genotyping
 prompting operator 11-5, 11-12
 specifying in .psd files 13-11
 STR genotyping
 .psd file 10-7
 plate setup template 9-4
 Plate Setup window 9-5
Sample Names tab
 sequencing 6-4, 6-5, 6-8
 STR genotyping 9-4, 9-5, 9-7
SCF files
 changing versions 5-10, 5-12
 MegaBACE.ini fields 5-4
 storage 15-4
sequencing
 administrator's task overview 1-3
 automatic base calling and file export 1-4, 5-7 to 5-13
 base callers, specifying 5-12 to 5-13
 Basecall.ini fields 5-11 to 5-12

data storage
 base called data 15-3 to 15-5
 raw data 15-1 to 15-2
file export options 5-10
file names and sample names 6-4, 6-5, 7-9, 15-1 to 15-2, 15-3 to 15-5
instrument control template
 creating 4-1 to 4-6
 overview 1-4
 setting default 5-5, 5-6
MegaBACE.ini fields 5-1 to 5-11
parameters, functional overview 1-7
plate definitions
 methods for creating 6-1 to 6-2
 workflow alternatives overview 1-5 to 1-6
plate setup data files (.psd)
 creating 7-4 to 7-6
 described 7-1
 field descriptions 6-4
 file name requirements 7-1 to 7-3
 workflow overview 1-5
plate setup template
 creating 6-3 to 6-8
 overview 1-4
 setting default 5-5 to 5-6, 6-9
Sequencing Diagnostics software annotation 5-12
service, serial numbers required xiv
SNP genotyping
 administrator's task overview 1-13
 data storage 15-2 to 15-3
 default application, specifying 11-5
 file names 12-4, 13-6, 15-2 to 15-3
 injections
 enforcing order of plates 11-3, 11-9
 maximum number 11-3, 11-7
 Sample Injection window 11-8
 specifying in master .psd file 13-9
 timeout period 11-3, 11-8
 instrument control template
 creating 4-1 to 4-6
 overview 1-14
 setting default 11-3, 11-7

-
- MegaBACE.ini fields 11-1 to 11-5
 - parameters, functional overview 1-18 to 1-19
 - plate attributes
 - .psd file 13-10 to 13-11
 - overview 1-18
 - prompting operator 11-4, 11-9 to 11-13
 - plate list, specifying in master .psd file 13-8
 - plate set definitions
 - methods for creating 12-1 to 12-2
 - workflow alternatives overview 1-15 to 1-16
 - plate set setup template
 - creating 12-4 to 12-8
 - overview 1-14
 - setting default 11-3, 11-6, 12-9
 - plate setup data files (.psd)
 - creating 13-4 to 13-5
 - described 13-1
 - sample names
 - .psd file 13-11
 - prompting operator 11-5, 11-12
 - SNP marker names
 - .psd file 13-10
 - prompting operator 11-4
 - SNP markers
 - .psd file, specifying in 13-6, 13-10
 - about software requirements 1-18
 - prompting operator for names 11-4, 11-9 to 11-12
 - software
 - Header Editor A-3
 - Instrument Control Manager 1-2
 - system components 1-1
 - Specify Plate ID window, illustrated 11-10
 - standard header format A-1
 - STR genotyping
 - administrator's task overview 1-8
 - data storage 15-1 to 15-2
 - default application, specifying 8-3
 - file names and sample names 9-4, 9-5, 10-7, 15-1 to 15-2
 - initialization files overview 1-9, 1-14
 - instrument control template
 - creating 4-1 to 4-6
 - overview 1-9
 - setting default 8-3, 8-4
 - MegaBACE.ini fields 8-1 to 8-3
 - parameters, functional overview 1-12
 - plate definitions
 - methods for creating 9-1 to 9-2
 - workflow alternatives overview 1-10 to 1-11
 - plate setup data files (.psd)
 - creating 10-4 to 10-6
 - described 10-1
 - field descriptions 9-4
 - file name requirements 10-1 to 10-3
 - plate setup template
 - creating 9-3 to 9-7
 - overview 1-9
 - setting default 8-3 to 8-4, 9-8
 - system initialization files. *see* Chemistry.ini file, Basecall.ini file, or MegaBACE.ini file
- T**
- tasks, administrator
 - sequencing 1-3
 - SNP genotyping 1-13
 - STR genotyping 1-8
 - templates, parameter. *see* instrument control templates (icp), plate setup templates (.tpl), or plate set setup templates
 - text files
 - creation of processed text files 5-5
 - creation of raw text files 5-5
- U**
- user documentation xiii
- W**
- workflows
 - sequencing, plate definition overview 1-5 to 1-6
 - SNP genotyping, plate set definition overview 1-15 to 1-16
 - STR genotyping, plate definition overview 1-10 to 1-11
-

Artisan Technology Group is an independent supplier of quality pre-owned equipment

Gold-standard solutions

Extend the life of your critical industrial, commercial, and military systems with our superior service and support.

We buy equipment

Planning to upgrade your current equipment? Have surplus equipment taking up shelf space? We'll give it a new home.

Learn more!

Visit us at artisantg.com for more info on price quotes, drivers, technical specifications, manuals, and documentation.

Artisan Scientific Corporation dba Artisan Technology Group is not an affiliate, representative, or authorized distributor for any manufacturer listed herein.

We're here to make your life easier. How can we help you today?

(217) 352-9330 | sales@artisantg.com | artisantg.com

210

Alamar Exhibit 1029
Alamar v. Olink



ARTISAN[®]
TECHNOLOGY GROUP