



Only Compatible with Solmetex® NXT Hg5® Amalgam Separator 1

Item# 26 98000

DRNA Generic Cartridge 2

Installation Instructions



You MUST replace your DRNA Generic Cartridge 2:

- When the sediment reaches the full line, or
- Every twelve months, whichever happens first.


Check the level of sedimentation on the DRNA Replacement Cartridge weekly

Included:


- Return Shipping Box/Foam
- (1) DRNA Replacement Cartridge
- (2) Yellow plugs to seal openings - Compatible with Solmetex® NXT Hg5® Amalgam Separator 1



PGR2023-00055 - Ex. 1005

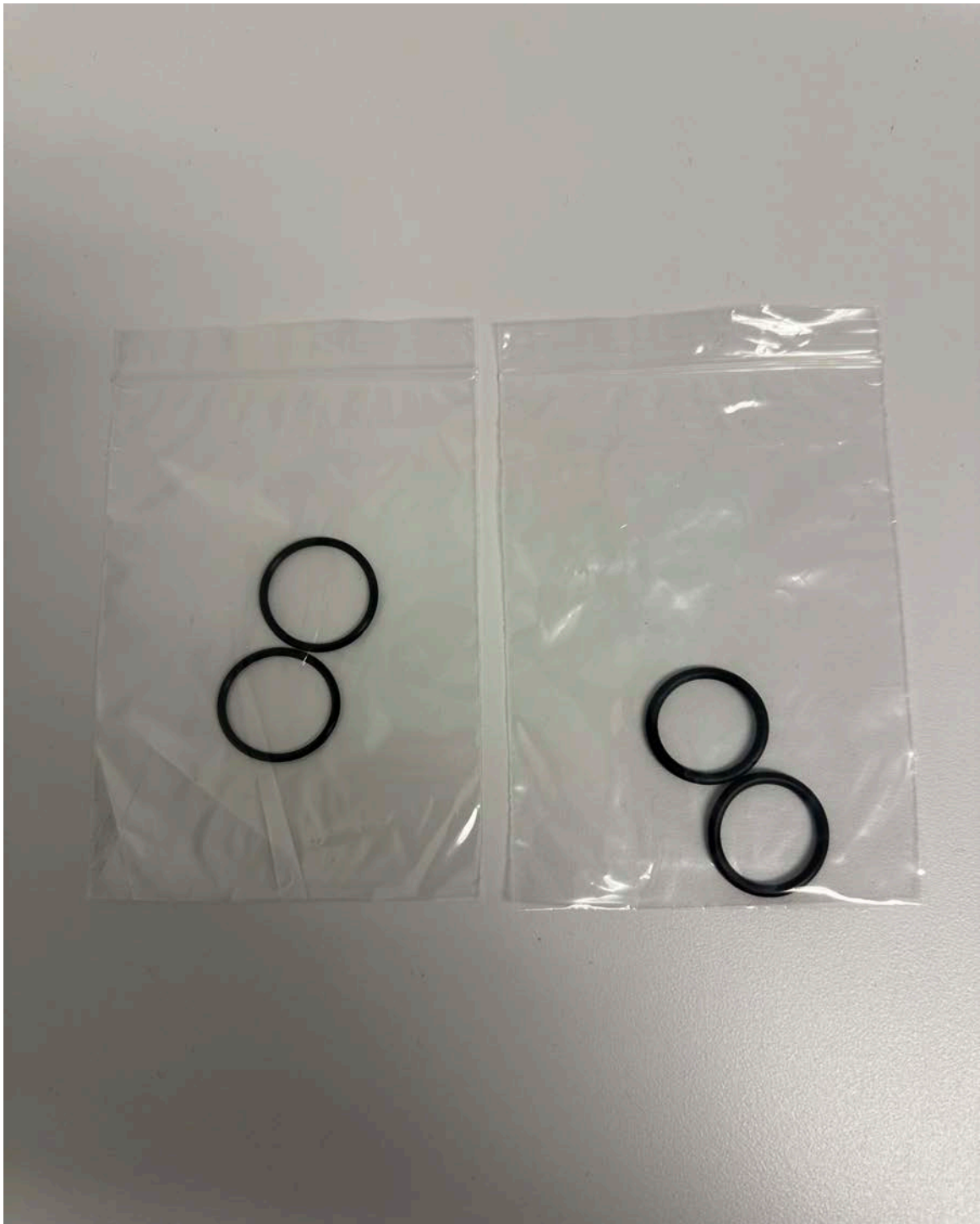

PGR2023-00055 - Ex. 1005

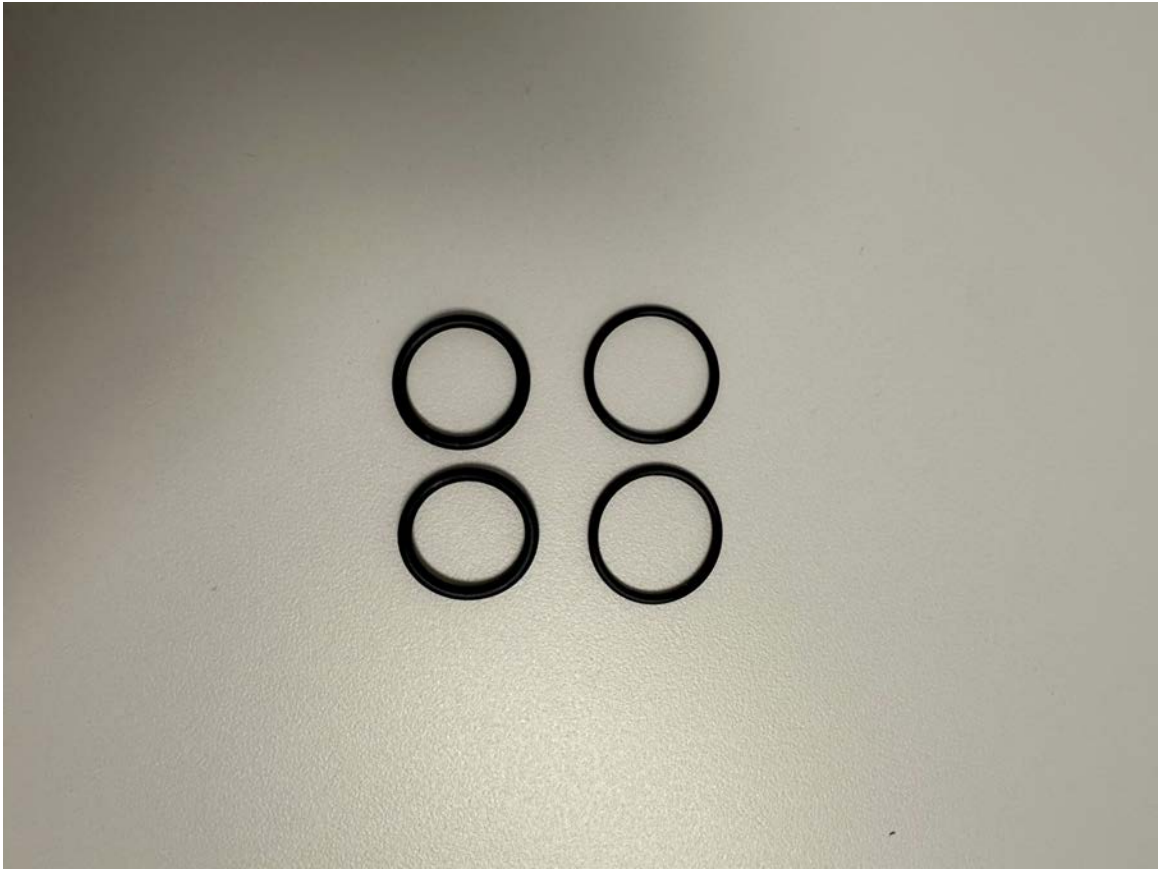

PGR2023-00055 - Ex. 1005



PGR2023-00055 - Ex. 1005



Ex. 1005



PGR2023-00055 - Ex. 1005



PGR2023-00055 - Ex. 1005



PGR2023-00055 - Ex. 1005



While the vacuum is turned ON and running:

Step 1: Remove full cartridge

- i. Open the zipper bag provided, remove and set aside (2) plugs to seal cartridge – depending on which cartridge you are returning, and slip bag underneath the installed cartridge.
- ii. While holding the installed cartridge with one hand, remove each of the two retaining pins using your other hand. Slowly pull cartridge down to release from the top compartment.
- iii. Locate the two openings on the top of the filled cartridge, and seal them with the (2) provided plugs. **Make sure that you securely push them all the way in, so that they cannot easily be removed.** You can use the palm of your hand to push down all the way.
- iv. With the filled cartridge inside the zipper bag, make sure to remove any surplus air and seal closed.

Step 2: Install DRNA replacement cartridge

- i. Wipe out port openings in top compartment with a paper towel or gauze. (This helps ensure a good seal)
- ii. Align the openings on the DRNA replacement cartridge with the openings underneath the upper compartment of your amalgam separator.
- iii. Push the DRNA replacement cartridge up and into the fittings on the top chamber. You may need to rock the cartridge gently until it securely connects to the top chamber and is seated.
- iv. Insert the retaining pins into their respective openings on each side of the DRNA replacement cartridge, until it is properly affixed to the top chamber.

Step 3: Fill out label on new DRNA Replacement Cartridge

- i. Write in the "Date Installed" date.
- ii. Write in the "Next Change Due" date (must be changed out 12 months from install date or once sediment reached the "Full Line" on the label, whichever occurs first).

Step 4: Check vacuum gauge

- i. Ensure that the reading on your vacuum gauge is at its normal level.

Step 5: Schedule Pickup

- i. Place used cartridge into the box that the new one came in – with top and bottom foam securing the cartridge in place within the box. Tape closed.
- ii. Contact DRNA at 800-360-1001 to request your prepaid return label:

Phone: 800-360-1001
Fax: 866-317-0853
Email: customerservice@drna.com

¹ Solmetex® and NXT Hg5® are registered trademarks of Solmetex LLC. DRNA is not affiliated with or sponsored or endorsed by Solmetex LLC.

Page 2 of 2

Dental Recycling North America, Inc. • 1270 Avenue of the Americas • New York, NY 10020 • Phone: 800-360-1001 • Fax: 866-317-0853 • Web: www.DRNA.com

The Leader in Dentalcare Waste Management!



Installation Instructions

DRNA Generic Cartridge 2

Item# 26 98000

Only Compatible with Solmetex®NXT Hg5® Amalgam Separator¹



You MUST replace your DRNA Generic Cartridge 2:

- **When the sediment reaches the full line, or**
- **Every twelve months, whichever happens first.**

Check the level of sedimentation on the DRNA Replacement Cartridge weekly



Included:

- Return Shipping Box/Foam
- (1) DRNA Replacement Cartridge
- (2) Yellow plugs to seal openings - *Compatible with Solmetex®NXT Hg5® Amalgam Separator¹*
- (2) Clear plugs to seal openings - *Compatible with DRNA Generic Cartridge 2*
- (1) Large Self-Sealing Zipper bag
- (2) Extra sets of O-rings - size medium and small. - *Usually not needed as cartridge comes with large O-rings already attached.*

Product Dimensions: 12.5" Height; 5.3" Diameter

Classification: Type 2, Sedimentation System

Max Fillable Volume: 1 Liter

Max Water Flow Rate: 1 Liter/minute

Line cleaners used must have a pH between 6 and 8

Prior to Starting Change Out Process:

- **Your vacuum system must be turned ON and running. You should keep it running for the duration of the change out process.**
- No patients are to be treated during the change out process.
- Make sure that there is no liquid in the upper compartment prior to performing change out. If there is liquid in upper compartment, postpone change out until the liquid has drained through to the filled lower collection filter previously installed.
- Wear medical gloves while performing change out.
- Check your vacuum gauge's normal level and take note of it.

¹ Solmetex® and NXT Hg5® are registered trademarks of Solmetex LLC. DRNA is not affiliated with or sponsored or endorsed by Solmetex LLC.

Page 1 of 2

Dental Recycling North America, Inc. • 1270 Avenue of the Americas • New York, NY 10020 • Phone: 800.360.1001 • Fax: 866.317.0853 • Web: www.DRNA.com


The Leader in Dentalcare Waste Management!!



PGR2023-00055 - Ex. 1005



PGR2023-00055 - Ex. 1005

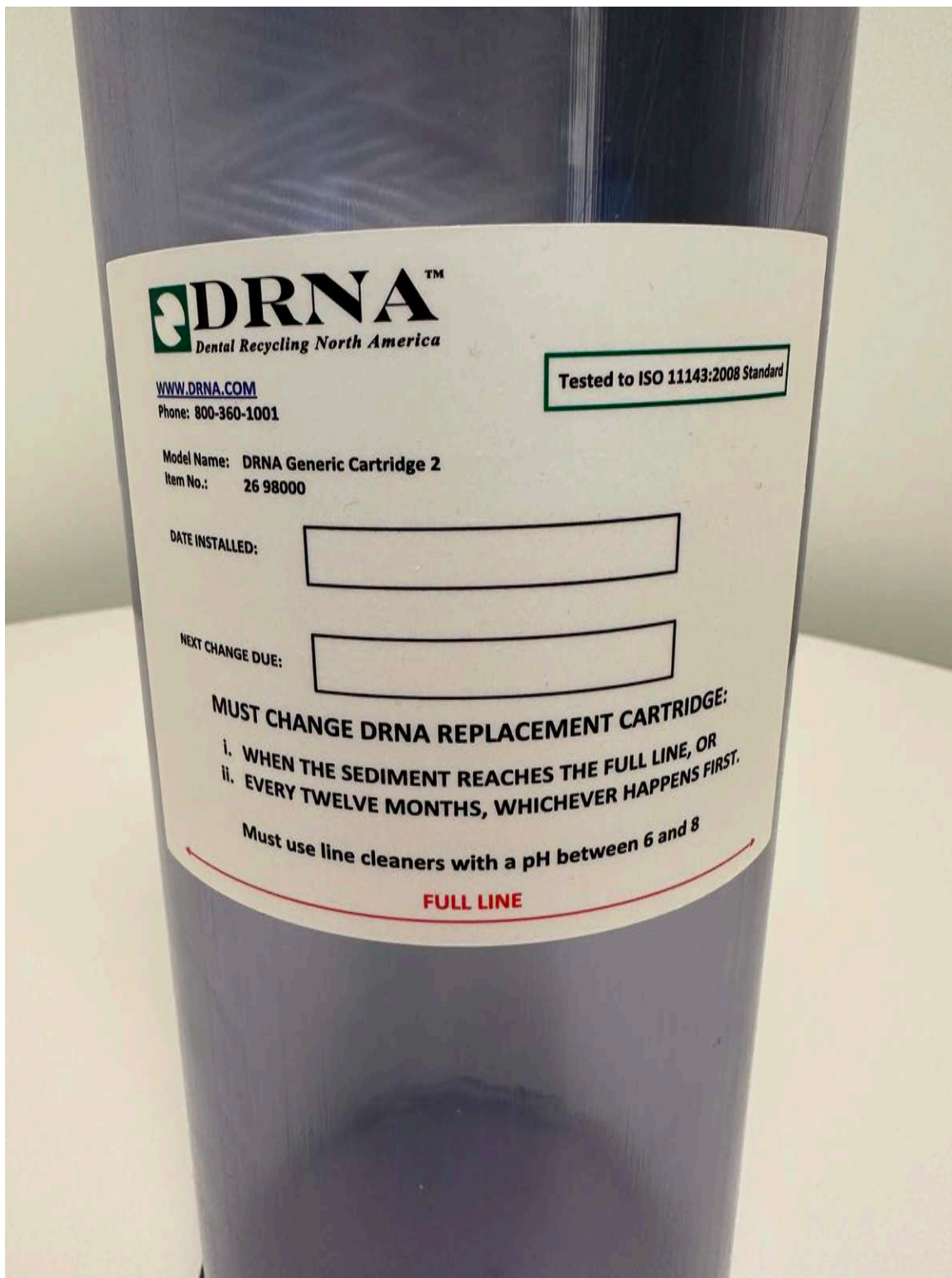

PGR2023-00055 - Ex. 1005
















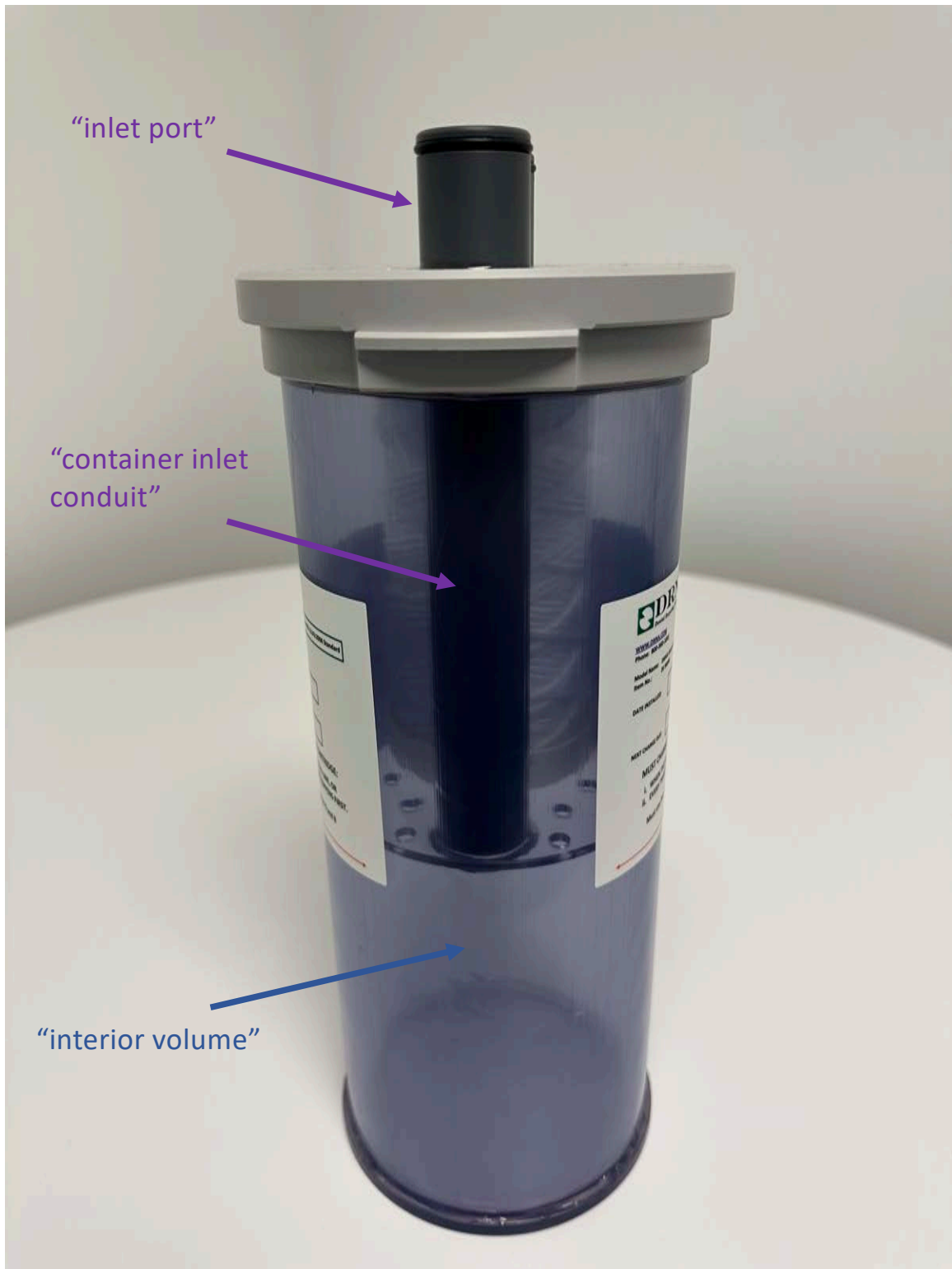





PGR2023-00055 - Ex. 1005




PGR2023-00055 - Ex. 1005



PGR2023-00055 - Ex. 1005



PGR2023-00055 - Ex. 1005

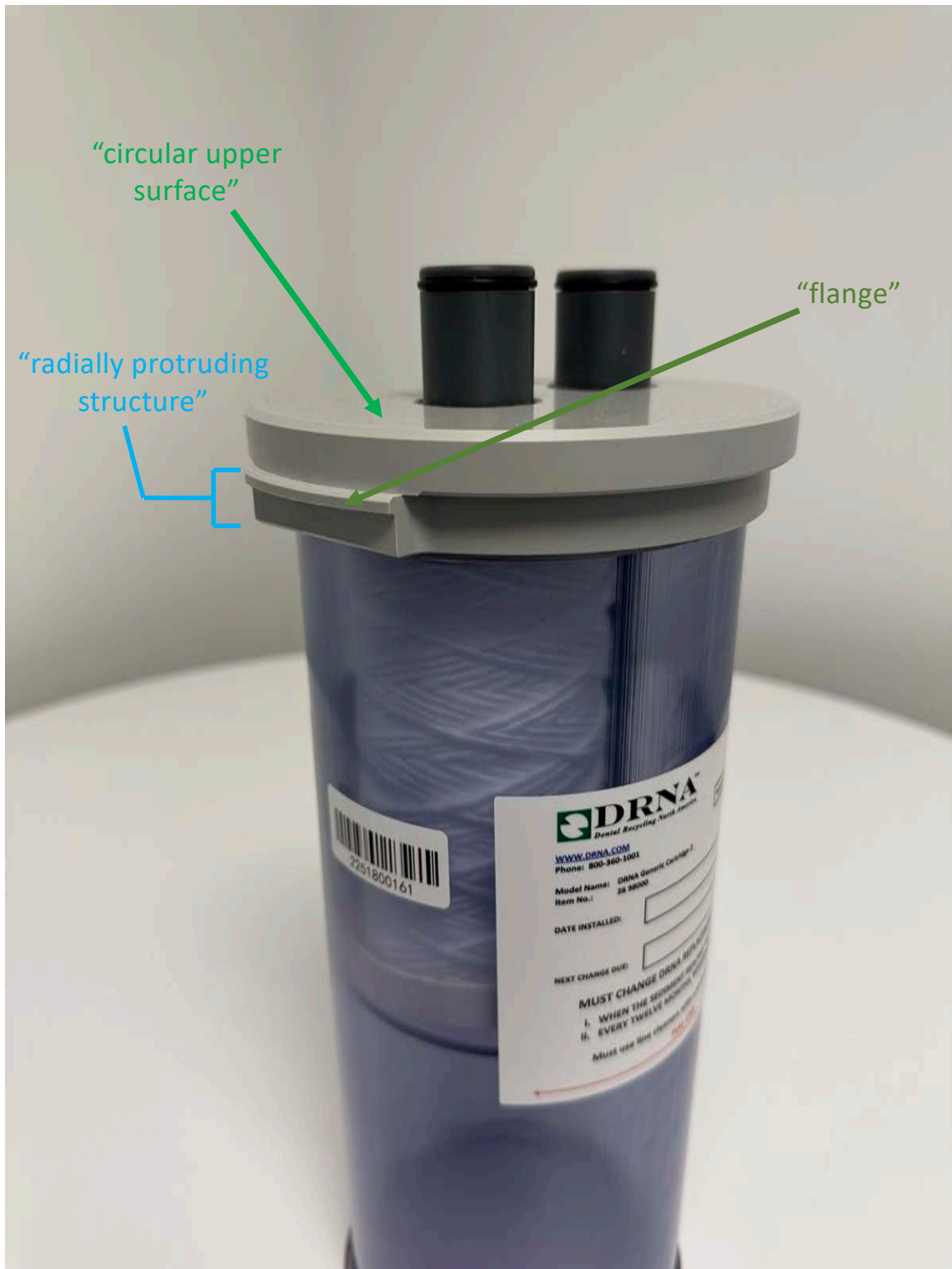

PGR2023-00055 - Ex. 1005





PGR2023-00055 - Ex. 1005



PGR2023-00055 - Ex. 1005



PGR2023-00055 - Ex. 1005



PGR2023-00055 - Ex. 1005