



Search

Create

9+



Downloads

1 / 1

Pin Guide & Fixation Base - CHROME Live Webinar Series (2/10)



ROE Dental Laboratory

1.2K views Apr 20, 2020

This webinar was recorded 3/25/2020, edited on 4/20. This webinar goes in-depth into all of the aspects of the Pin Guide and the Fixation Base. It also includes Q&A from webinar attendees. For more information on CHROME please visit our website - https://roedentallab.com/CHROME_YT Subscribe to ROE Dental Laboratory on YouTube: <https://bit.ly/38bjrPy> Like ROE on Facebook: <https://www.facebook.com/roedental> Follow ROE on LinkedIn: <https://www.linkedin.com/in/roe-dental-laboratory/> Follow ROE on Instagram: <https://www.instagram.com/roedentallab/> Follow ROE on Twitter: <https://twitter.com/roedentallab> To receive news, product updates, event notices etc: <https://www.roedentallab.com/news> For more information, stay tuned to: <https://www.roedentallab.com/>

Show less



CHROME Live Webinar Series

The Pin Guide and Fixation Base



Pin Guide – Inspection

Copyright ROE Dental Laboratory

Dentate:


Fits the model, will fit the teeth

Edentulous:

Fits the patient like the denture – it is a duplicate of the denture

Loops should Passively fit the CHROME Locs

DO NOT seat PG to the Fixation Base crooked – the loops will break.

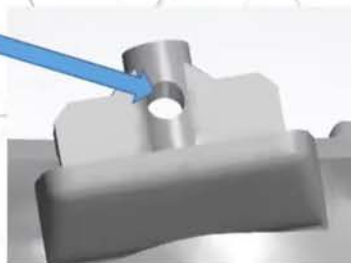

CHROME™

FULL ARCH GUIDED STABILITY

Edentulous Troubleshooting

Copyright ROE Dental Laboratory



CHROME™
FULL ARCH GUIDED STABILITY



Copyright ROE Dental Laboratory

Pin Guide Single Edentulous – Guide will only fit the mouth

- Will fit just like the patient's denture. If it rocks, adjust like a denture. Will meet the opposing teeth. Patient in occlusion while drilling facial pins.

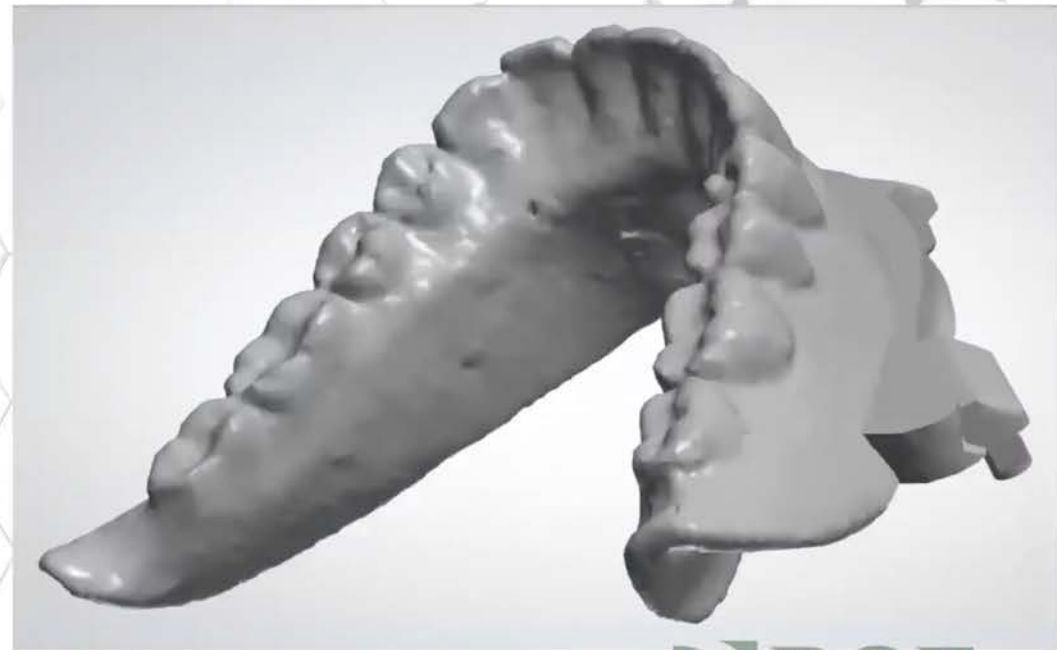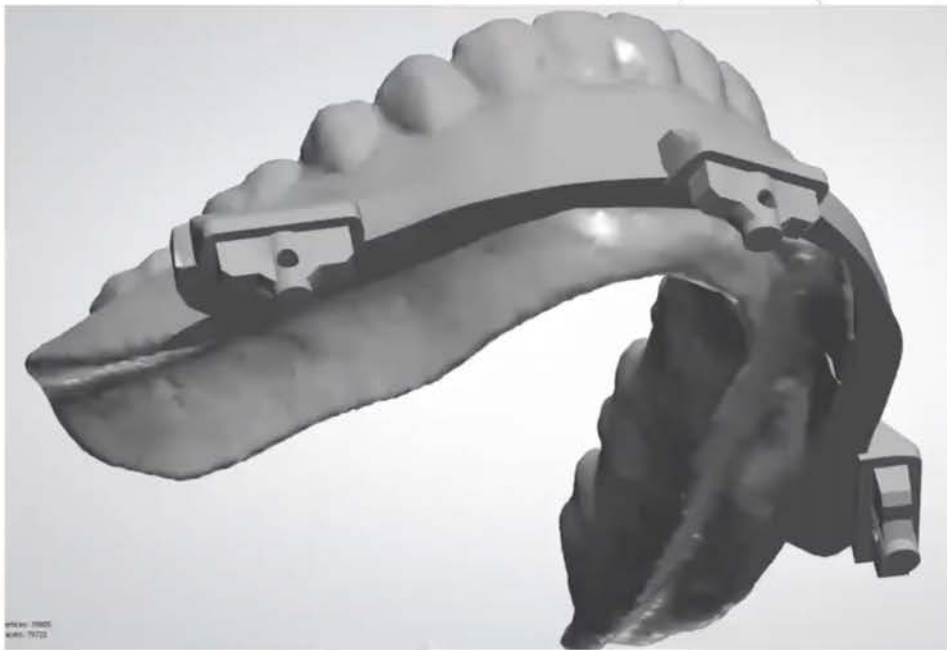

CHROME
FULL ARCH GUIDED STABILITY



ROE
Dental Laboratory

Single Edentulous integrated bite



Copyright ROE Dental Laboratory

Copyright ROE Dental Laboratory

Pin Guide Double Edentulous – U/L do not fit together

- Maxillary Pin Guide fits against existing denture
- Mandibular Pin Guide fits against new maxillary prosthetic
- Will deviate from this if the set-up is different from the dentures

