

Final Report

Antibody Full Length Sequencing of Hybridoma 5G1.1 T175

Order Number: U856UEK140-1

November 15, 2019

Prepared by

Jinquan Chen

Jinquan Chen, Production Supervisor

This study was conducted by strictly following the procedures described in this report. All data presented are authentic and accurate to the best of our knowledge. Information and materials provided by this investigation are for research use only.

Recommended Analysis Tools

We recommend the following free online tools for sequence analysis of immunoglobulin variable regions:

- i. NCBI Nucleotide BLAST
- ii. IMGT/V Quest program
- iii. NCBI IgBLAST

Antibody sequences of 5G1.1 T175 were listed as below:

Heavy chain: DNA sequence (1416 bp)

Signal sequence-**FR1-CDR1-FR2-CDR2-FR3-CDR3-FR4-Constant region-Stop codon**

ATGGAATGGACCTGGGTCTTTCTCTTCCTCCTGTCAGTAACTGCAGGTGTCCACTCC**CAGGTT**CAGC
TGCAGCAGTCTGGAGCTGAGCTGATGAAGCCTGGGGCCTCAGTGAAGATGTCCTGCAAGGCTACTGG
CTACATATTCAGT**AACTACTGGATACAGTGGATAAAGCAGAGGCCTGGACATGGCCTTGAGTGGATT**
GGTGAGATTTT**ACCTGGAAGTGGTTCTACTGAGTACACTGAGA**ACTT**CAAGGACAAGGCCGCATTCA**
CTGCAGATACATCCTCCAACACAGCCTACATGCAACTCAGCAGCCTGACATCAGAGGACTCTGCCGT
CTATTACTGTGCAAGAT**TATTTCTTCGGTAGTAGCCCCAACTGGTACTTCGATGTC**TGGGGCGCAGGG
ACCACGGTCACCGTCTCCT**CAGCCAAAACAACAGCCCCATCGGTCTATCCACTGGCCCCCTGTGTGTG**
GAGATACA**ACTGGCTCCTCGGTGACTCTAGGATGCCTGGTCAAGGTTATTTCCCTGAGCCAGTGAC**
CTTGACCTGGA**ACTCTGGATCCCTGTCCAGTGGTGTGCACACCTTCCCAGCTGTCCTGCAGTCTGAC**
CTCTACAC**CTCAGCAGCTCAGTGACTGTAACCTCGAGCACCTGGCCCAGCCAGTCCATCACCTGCA**
ATGTGGCC**ACCCGGCAAGCAGCACCAAGGTGGACAAGAAAATTGAGCCCAGAGGGCCCACAATCAA**
GCCCTGTCCTCCATGCAAATGCCCAGCACCTAACCTCTTGGGTGGACCATCCGTCTTCATCTTCCCT
CCAAAGATCAAGGATGTACTCATGATCTCCCTGAGCCCCATAGTCACATGTGTGGTGGTGGATGTGA
GCGAGGATGACCCAGATGTCCAGATCAGCTGGTTTGTGAACAACGTGGAAGTACACACAGCTCAGAC
ACAAACCCATAGAGAGGATTACAACAGTACTCTCCGGGTGGTCAGTGCCCTCCCCATCCAGCACCAG
GACTGGATGAGTGGCAAGGAGTTCAAATGCAAGGTCAACAACAAAGACCTCCCAGCGCCCATCGAGA
GAACCATCTCAAAACCCAAAGGGTCAGTAAGAGCTCCACAGGTATATGTCTTGCCCTCCACCAGAAGA
AGAGATGACTAAGAAACAGGTC**ACTCTGACCTGCATGGTCACAGACTTCATGCCTGAAGACATTTAC**
GTGGAGTGGACCAACAACGGGAAAACAGAGCTAAACTACAAGA**AACTGAACCAGTCCTGGACTCTG**
ATGGTTCTTACTT**CATGTACAGCAAGCTGAGAGTGGA**AAAGAAGAACTGGGTGGAAAGAAATAGCTA
CTCCTGTT**CAGTGGTCCACGAGGGTCTGCACAATCACCACACGACTAAGAGCTTCTCCCGGACTCCG**
GGTAAATGA

Heavy chain: Amino acid sequence (471 aa)

Signal peptide-**FR1-CDR1-FR2-CDR2-FR3-CDR3-FR4-Constant region-Stop codon**

MEWTWVFLFLLSVTAGVHS**QVQLQQSGAELMKPGASVKMSCKATGYIFS****NYWIQ**WIKQRP**GHGLEWI**

GEILPGSGSTEYTENFKDKAAFTADTSSNTAYMQLSSLTSEDSAVYYCARYFFGSSPNWYFDVWGAG
TTVTVSSAKTTAPSVYPLAPVCGD TTGSSVTLGCLVKGYFPEPVTLTWNSGSLSSGVHTFPAVLQSD
LYTLSSSVTVTSSTWPSQSITCNVAHPASSTKVDKKIEPRGPTIKPCPPCKCPAPNLLGGPSVFIFP
PKIKDVLMLISLSPIVTCVVVDVSEDDPDVQISWVFVNNVEVHTAQQTQTHREDYNSTLRVVSALPIQHQ
DWMSGKEFKCKVNNKDLPAPIERTISKPKGVRAPQVYVLPPEEEMTKKQVTLTCMVTD FMPEDIY
VEWTNNGKTELNYKNTEPVLDSG SYFMYSKLRVEKKNWVERNSYSCSVVHEGLHNHHTTKSFSRTP
GK-

Light chain: DNA sequence (702 bp)

Signal sequence-FR1-CDR1-FR2-CDR2-FR3-CDR3-FR4-Constant region-Stop codon

ATGGGTGTACCCACTCAGCTCCTGTTGCTGTGGCTTACAGTCGTAATTGTCAGATGTGACATCCAGA
TGACTCAGTCTCCAGCTTCACTGTCTGCATCTGTGGGAGAACTGTCACCATCACATGTGGAGCAAG
TGAGAATATTTACGGTGCTTTAAATTGGTATCAGCGGAAACAGGGAAAATCTCCTCAGCTCCTGATC
TATGGTGCAACCAACTTGGCAGATGGCATGTCATCGAGGTTCACTGGCAGTGGATCTGGTAGACAGT
ATTATCTCAAGATCAGTAGCCTGCATCCTGACGATGTTGCAACGTATTACTGTCAAATGTGTTAAA
TACTCCTCTCACGTTTCGGTGCTGGGACCAAGTTGGAGCTGAAACGGGCTGATGCTGCACCAACTGTA
TCCATCTTCCCACCATCCAGTGAGCAGTTAACATCTGGAGGTGCCTCAGTCGTGTGCTTCTTGAACA
ACTTCTACCCCAAAGACATCAATGTCAAGTGGAAGATTGATGGCAGTGAACGACAAAATGGCGTCCT
GAACAGTTGGACTGATCAGGACAGCAAAGACAGCACCTACAGCATGAGCAGCACCTCACGTTGACC
AAGGACGAGTATGAACGACATAACAGCTATACCTGTGAGGCCACTCACAAGACATCAACTTCACCCA
TTGTCAAGAGCTTCAACAGGAATGAGTGT TAG


Light chain: Amino acid sequence (233 aa)

Signal peptide-FR1-CDR1-FR2-CDR2-FR3-CDR3-FR4-Constant region-Stop codon

MGVPTQLLLLWLTVVIVRC DIQMTQSPASLSASVGETVTITCGASENIYGALN WYQRKQKSPQLLI
YGATNLADGMSSRFSGSGRQYYLKISSLHPDDVATYYCQNVLNTPLTFGAGTKLELKRADAAPTV
SIFPPSSEQLTSGGASVVCFLNNFY PKDINVKWKIDGSERQNGVLNSWTDQDSK DSTYSMSSTLTLT
KDEYERHNSYTCEATHKTSTSPIVKSFNRNEC-

Project Specifications

Sequence	Clones sequenced	Clones with >99% sequence identity
V _H	5	5
C _H	5	5
V _L	5	5
C _L	5	5

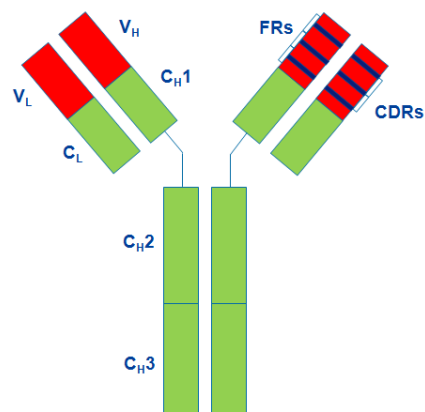

IMGT Analysis of V(D)J Junctions

Sequence	V-GENE and allele	Functionality	V-REGION identity % (nt)	J-GENE and allele	D-GENE and allele	AA Junction	Junction frame
V _H	Musmus IGHV1-9*01 F	productive	93.06% (268/288 nt)	Musmus IGHJ1*01 F	Musmus IGHD1-1*01 F	CARYFF GSSPNW YFDVW	in-frame
V _L	Musmus IGKV12-89*01 F	productive	99.28% (277/279 nt)	Musmus IGKJ5*01 F		CQNVLN TPLTF	in-frame

Isotype of antibody

Isotype is mouse IgG2a/kappa, as analyzed by the sequences of constant region.

Materials

Hybridoma cells provided by the client;

TRIzol® Reagent (Ambion, Cat. No. : 15596-026);

PrimeScript™ 1st Strand cDNA Synthesis Kit (Takara, Cat. No. : 6110A).

Methods

Total RNA was isolated from the hybridoma cells following the technical manual of TRIzol® Reagent. Total RNA was then reverse-transcribed into cDNA using either isotype-specific anti-sense primers or universal primers following the technical manual of PrimeScript™ 1st Strand cDNA Synthesis Kit. Antibody fragments of heavy chain and light chain were amplified according to the standard operating procedure (SOP) of rapid amplification of cDNA ends (RACE) of GenScript. Amplified antibody fragments were cloned into a standard cloning vector separately. Colony PCR was performed to screen for clones with inserts of correct sizes. The consensus sequence was provided.