



IEC 60601-2-2

Edition 5.0 2009-02

INTERNATIONAL STANDARD

NORME INTERNATIONALE

**Medical electrical equipment –
Part 2-2: Particular requirements for the basic safety and essential performance
of high frequency surgical equipment and high frequency surgical accessories**

**Appareils électromédicaux –
Partie 2-2: Exigences particulières pour la sécurité de base et les performances
essentiels des appareils d'électrochirurgie à courant haute fréquence et des
accessoires d'électrochirurgie à courant haute fréquence**

IEC 60601-2-2:2009

Customer: Bradley Hyde - No. of User(s): 1 - Company: Latham and Watkins
Order No.: WS-2023-008299 - IMPORTANT: This file is copyright of IEC, Geneva, Switzerland. All rights reserved.
This file is subject to a licence agreement. Enquiries to Email: sales@iec.ch - Tel.: +41 22 919 02 11

SERENDIA-00084149

CX-2848.1

Serendia, LLC
Exhibit 2077
EndyMed v. Serendia
IPR2024-00842
Ex. 2077-001



THIS PUBLICATION IS COPYRIGHT PROTECTED

Copyright © 2009 IEC, Geneva, Switzerland

All rights reserved. Unless otherwise specified, no part of this publication may be reproduced or utilized in any form or by any means, electronic or mechanical, including photocopying and microfilm, without permission in writing from either IEC or IEC's member National Committee in the country of the requester.

If you have any questions about IEC copyright or have an enquiry about obtaining additional rights to this publication, please contact the address below or your local IEC member National Committee for further information.

Droits de reproduction réservés. Sauf indication contraire, aucune partie de cette publication ne peut être reproduite ni utilisée sous quelque forme que ce soit et par aucun procédé, électronique ou mécanique, y compris la photocopie et les microfilms, sans l'accord écrit de la CEI ou du Comité national de la CEI du pays du demandeur.

Si vous avez des questions sur le copyright de la CEI ou si vous désirez obtenir des droits supplémentaires sur cette publication, utilisez les coordonnées ci-après ou contactez le Comité national de la CEI de votre pays de résidence.

IEC Central Office
3, rue de Varembe
CH-1211 Geneva 20
Switzerland
Email: inmail@iec.ch
Web: www.iec.ch

About the IEC

The International Electrotechnical Commission (IEC) is the leading global organization that prepares and publishes International Standards for all electrical, electronic and related technologies.

About IEC publications

The technical content of IEC publications is kept under constant review by the IEC. Please make sure that you have the latest edition, a corrigenda or an amendment might have been published.

■ Catalogue of IEC publications: www.iec.ch/searchpub

The IEC on-line Catalogue enables you to search by a variety of criteria (reference number, text, technical committee,...). It also gives information on projects, withdrawn and replaced publications.

■ IEC Just Published: www.iec.ch/online_news/justpub

Stay up to date on all new IEC publications. Just Published details twice a month all new publications released. Available on-line and also by email.

■ Electropedia: www.electropedia.org

The world's leading online dictionary of electronic and electrical terms containing more than 20 000 terms and definitions in English and French, with equivalent terms in additional languages. Also known as the International Electrotechnical Vocabulary online.

■ Customer Service Centre: www.iec.ch/webstore/custserv

If you wish to give us your feedback on this publication or need further assistance, please visit the Customer Service Centre FAQ or contact us:

Email: csc@iec.ch
Tel.: +41 22 919 02 11
Fax: +41 22 919 03 00

A propos de la CEI

La Commission Electrotechnique Internationale (CEI) est la première organisation mondiale qui élabore et publie des normes internationales pour tout ce qui a trait à l'électricité, à l'électronique et aux technologies apparentées.

A propos des publications CEI

Le contenu technique des publications de la CEI est constamment revu. Veuillez vous assurer que vous possédez l'édition la plus récente, un corrigendum ou amendement peut avoir été publié.

■ Catalogue des publications de la CEI: www.iec.ch/searchpub/cur_fut-f.htm

Le Catalogue en-ligne de la CEI vous permet d'effectuer des recherches en utilisant différents critères (numéro de référence, texte, comité d'études,...). Il donne aussi des informations sur les projets et les publications retirées ou remplacées.

■ Just Published CEI: www.iec.ch/online_news/justpub

Restez informé sur les nouvelles publications de la CEI. Just Published détaille deux fois par mois les nouvelles publications parues. Disponible en-ligne et aussi par email.

■ Electropedia: www.electropedia.org

Le premier dictionnaire en ligne au monde de termes électroniques et électriques. Il contient plus de 20 000 termes et définitions en anglais et en français, ainsi que les termes équivalents dans les langues additionnelles. Egalement appelé Vocabulaire Electrotechnique International en ligne.

■ Service Clients: www.iec.ch/webstore/custserv/custserv_entry-f.htm

Si vous désirez nous donner des commentaires sur cette publication ou si vous avez des questions, visitez le FAQ du Service clients ou contactez-nous:

Email: csc@iec.ch
Tél.: +41 22 919 02 11
Fax: +41 22 919 03 00

Customer: Bradley Hyde- No. of User(s): 1 - Company: Latham and Watkins

Order No.: WS-2023-008299 - IMPORTANT: This file is copyright of IEC, Geneva, Switzerland. All rights reserved.

This file is subject to a licence agreement. Enquiries to Email: sales@iec.ch - Tel.: +41 22 919 02 11

SERENDIA-00084150

Ex. 2077-002

CX-2848.2

INTERNATIONAL STANDARD

NORME INTERNATIONALE

**Medical electrical equipment –
Part 2-2: Particular requirements for the basic safety and essential performance
of high frequency surgical equipment and high frequency surgical accessories**

**Appareils électromédicaux –
Partie 2-2: Exigences particulières pour la sécurité de base et les performances
essentielle des appareils d'électrochirurgie à courant haute fréquence et des
accessoires d'électrochirurgie à courant haute fréquence**

INTERNATIONAL
ELECTROTECHNICAL
COMMISSION

COMMISSION
ELECTROTECHNIQUE
INTERNATIONALE

PRICE CODE
CODE PRIX **XC**

ICS 11.040.30

ISBN 978-2-88910-213-6

Customer: Bradley Hyde, No. of Users: 1, Company: Latham and Watkins
is a registered trademark of the International Electrotechnical Commission.
Order No. 4430023-000289 - IEC 60601-2-2:2009-02
This file is subject to a licence agreement. Enquiries to Email: sales@iec.ch - Tel.: +41 22 919 02 11

SERENDIA-00084151

Ex. 2077-003

CX-2848.3

CONTENTS

FOREWORD.....	4
INTRODUCTION.....	7
201.1 Scope, object and related standards	8
201.2 Normative references	9
201.3 Terms and definitions.....	10
201.4 General requirements.....	13
201.5 General requirements for testing of ME EQUIPMENT.....	14
201.6 Classification of ME EQUIPMENT and ME SYSTEMS	14
201.7 ME EQUIPMENT identification, marking and documents.....	14
201.8 Protection against electrical HAZARDS from ME EQUIPMENT.....	19
201.9 Protection against mechanical hazards of ME EQUIPMENT and ME SYSTEMS.....	33
201.10 Protection against unwanted and excessive radiation HAZARDS	33
201.11 Protection against excessive temperatures and other HAZARDS.....	34
201.12 Accuracy of controls and instruments and protection against hazardous outputs.....	35
201.13 HAZARDOUS SITUATIONS and fault conditions.....	41
201.14 PROGRAMMABLE ELECTRICAL MEDICAL SYSTEMS (PEMS).....	42
201.15 Construction of ME EQUIPMENT	42
201.16 ME SYSTEMS	46
201.17 Electromagnetic compatibility of ME EQUIPMENT and ME SYSTEMS	46
202 * Electromagnetic compatibility – Requirements and tests.....	46
208 General requirements, tests and guidance for alarm systems in medical electrical equipment and medical electrical systems.....	47
Annexes	47
Annex AA (informative) Particular guidance and rationale.....	48
Annex BB (informative) ELECTROMAGNETIC DISTURBANCES created by HF SURGICAL EQUIPMENT.....	70
Bibliography.....	79
Figure 201.101 – Symbol used with an EARTH REFERENCED PATIENT CIRCUIT.....	14
Figure 201.102 – Symbol used with a HF ISOLATED PATIENT CIRCUIT	15
Figure 201.103 – Circuit suitable for testing compliance to 201.8.4.101	20
Figure 201.104 – Measurement of HF LEAKAGE CURRENT with NEUTRAL ELECTRODE referenced to earth and load between electrodes.....	23
Figure 201.105 – Measurement of HF LEAKAGE CURRENT with NEUTRAL ELECTRODE referenced to earth and load from ACTIVE ELECTRODE to earth	24
Figure 201.106 – Measurement of HF LEAKAGE CURRENT with NEUTRAL ELECTRODE isolated from earth at HIGH FREQUENCY.....	25
Figure 201.107 – Measurement of HF LEAKAGE CURRENT from a BIPOLAR ELECTRODE	26
Figure 201.108 – Test apparatus for anchorages of cords of ACTIVE ACCESSORY.....	32
Figure 201.109 – Measurement of output power – MONOPOLAR output.....	37
Figure 201.110 – Measurement of output power – BIPOLAR output.....	38

Customer: Bradley Hyde- No. of User(s): 1 - Company: Latham and Watkins

Order No.: WS-2023-008299 - IMPORTANT: This file is copyright of IEC, Geneva, Switzerland. All rights reserved.

This file is subject to a licence agreement. Enquiries to Email: sales@iec.ch - Tel.: +41 22 919 02 11

Figure 201.111 – Method of testing feedback from one active output to another in simultaneous activation.....	41
Figure AA.1 – Example of various parts of an HF surgical system.....	49
Figure AA.2 – CREST FACTOR vs. peak voltage	53
Figure AA.3 – Example of PATIENT circuit with NEUTRAL ELECTRODE referenced to earth at operating frequencies	57
Figure BB.1 – E-FIELD EMISSIONS test setup.....	73
Figure BB.2 – H-FIELD EMISSIONS test setup	74
Figure BB.3 – Conducted EMISSIONS test setup	75
Figure BB.4 – Unit ad hoc test	77
Figure BB.5 – Power cord ad hoc test.....	78
Figure BB.6 – ACCESSORY cord ad hoc test	78
Table 201.101 – Colours of indicator lights and their meaning for HF SURGICAL EQUIPMENT	15
Table 201.102 – Maximum output powers in SINGLE FAULT CONDITIONS	40
Table 201.103 – Test currents by weight range.....	44
Table AA.1 – Summary of measured current and durations for 25 TUR procedures.....	65
Table AA.2 – Summary of measured currents and durations for general surgical procedures	66
Table BB.1 – Worst case emissions of spark gap type HF SURGICAL EQUIPMENT	76
Table BB.2 – Worst case emissions of non-spark gap (modern) HF SURGICAL EQUIPMENT	76

INTERNATIONAL ELECTROTECHNICAL COMMISSION

MEDICAL ELECTRICAL EQUIPMENT –

Part 2-2: Particular requirements for the basic safety and essential performance of high frequency surgical equipment and high frequency surgical accessories

FOREWORD

- 1) The International Electrotechnical Commission (IEC) is a worldwide organization for standardization comprising all national electrotechnical committees (IEC National Committees). The object of IEC is to promote international co-operation on all questions concerning standardization in the electrical and electronic fields. To this end and in addition to other activities, IEC publishes International Standards, Technical Specifications, Technical Reports, Publicly Available Specifications (PAS) and Guides (hereafter referred to as "IEC Publication(s)"). Their preparation is entrusted to technical committees; any IEC National Committee interested in the subject dealt with may participate in this preparatory work. International, governmental and non-governmental organizations liaising with the IEC also participate in this preparation. IEC collaborates closely with the International Organization for Standardization (ISO) in accordance with conditions determined by agreement between the two organizations.
- 2) The formal decisions or agreements of IEC on technical matters express, as nearly as possible, an international consensus of opinion on the relevant subjects since each technical committee has representation from all interested IEC National Committees.
- 3) IEC Publications have the form of recommendations for international use and are accepted by IEC National Committees in that sense. While all reasonable efforts are made to ensure that the technical content of IEC Publications is accurate, IEC cannot be held responsible for the way in which they are used or for any misinterpretation by any end user.
- 4) In order to promote international uniformity, IEC National Committees undertake to apply IEC Publications transparently to the maximum extent possible in their national and regional publications. Any divergence between any IEC Publication and the corresponding national or regional publication shall be clearly indicated in the latter.
- 5) IEC provides no marking procedure to indicate its approval and cannot be rendered responsible for any equipment declared to be in conformity with an IEC Publication.
- 6) All users should ensure that they have the latest edition of this publication.
- 7) No liability shall attach to IEC or its directors, employees, servants or agents including individual experts and members of its technical committees and IEC National Committees for any personal injury, property damage or other damage of any nature whatsoever, whether direct or indirect, or for costs (including legal fees) and expenses arising out of the publication, use of, or reliance upon, this IEC Publication or any other IEC Publications.
- 8) Attention is drawn to the Normative references cited in this publication. Use of the referenced publications is indispensable for the correct application of this publication.
- 9) Attention is drawn to the possibility that some of the elements of this IEC Publication may be the subject of patent rights. IEC shall not be held responsible for identifying any or all such patent rights.

International standard IEC 60601-2-2 has been prepared by IEC subcommittee 62D: Electromedical equipment, of IEC technical committee 62: Electrical equipment in medical practice.

This fifth edition cancels and replaces the fourth edition published in 2006. This edition constitutes a technical revision. Revisions in this edition include new language for preconditioning accessories prior to insulation testing, refining the requirements for electromagnetic compatibility testing and correcting some of the equations used in deriving the thermal test for NEUTRAL ELECTRODES.

Customer: Bradley Hyde- No. of User(s): 1 - Company: Latham and Watkins

Order No.: WS-2023-008299 - IMPORTANT: This file is copyright of IEC, Geneva, Switzerland. All rights reserved.

This file is subject to a licence agreement. Enquiries to Email: sales@iec.ch - Tel.: +41 22 919 02 11

SERENDIA-00084154

Ex. 2077-006

CX-2848.6

The text of this particular standard is based on the following documents:

FDIS	Report on voting
62D/726/FDIS	62D/755/RVD

Full information on the voting for the approval of this particular standard can be found in the report on voting indicated in the above table.

This publication has been drafted in accordance with the ISO/IEC Directives, Part 2.

In this standard, the following print types are used:

- Requirements and definitions: roman type.
- *Test specifications: italic type.*
- Informative material appearing outside of tables, such as notes, examples and references, in smaller type. Normative text of tables is also in a smaller type.
- TERMS DEFINED IN CLAUSE 3 OF THE GENERAL STANDARD, IN THIS PARTICULAR STANDARD OR AS NOTED: SMALL CAPITALS.

In referring to the structure of this standard, the term

- “clause” means one of the seventeen numbered divisions within the table of contents, inclusive of all subdivisions (e.g. Clause 7 includes subclauses 7.1, 7.2, etc.);
- “subclause” means a numbered subdivision of a clause (e.g. 7.1, 7.2 and 7.2.1 are all subclauses of Clause 7).

References to clauses within this standard are preceded by the term “Clause” followed by the clause number. References to subclauses within this standard are by number only.

In this standard, the conjunctive “or” is used as an “inclusive or” so a statement is true if any combination of the conditions is true.

The verbal forms used in this standard conform to usage described in Annex H of the ISO/IEC Directives, Part 2. For the purposes of this standard, the auxiliary verb:

- “shall” means that compliance with a requirement or a test is mandatory for compliance with this standard;
- “should” means that compliance with a requirement or a test is recommended but is not mandatory for compliance with this standard;
- “may” is used to describe a permissible way to achieve compliance with a requirement or test.

An asterisk (*) as the first character of a title or at the beginning of a paragraph or table title indicates that there is guidance or rationale related to that item in Annex AA.

A list of all parts of the IEC 60601 series, published under the general title *Medical electrical equipment*, can be found on the IEC website.

The committee has decided that the contents of this particular standard will remain unchanged until the maintenance result date indicated on the IEC web site under "<http://webstore.iec.ch>" in the data related to the specific publication. At this date, the publication will be

- reconfirmed;
- withdrawn;
- replaced by a revised edition, or
- amended

The contents of the corrigendum of February 2014 have been included in this copy.

Withdrawn

INTRODUCTION

The minimum safety requirements specified in this particular standard are considered to provide for a practical degree of safety in the operation of high frequency surgical equipment.

This particular standard amends and supplements IEC 60601-1 (third edition, 2005): *Medical electrical equipment – Part 1: General requirements for basic safety and essential performance*, hereinafter referred to as the general standard (see 201.1.4).

The requirements are followed by specifications for the relevant tests.

A "Particular guidance and rationale" section giving some explanatory notes, where appropriate, about the more important requirements is included in annex AA.

Clauses or subclauses for which there are explanatory notes in annex AA are marked with an asterisk (*).

It is considered that a knowledge of the reasons for these requirements will not only facilitate the proper application of the standard but will, in due course, expedite any revision necessitated by changes in clinical practice or as a result of developments in technology. However, this annex does not form part of the requirements of this standard.

MEDICAL ELECTRICAL EQUIPMENT –

Part 2-2: Particular requirements for the basic safety and essential performance of high frequency surgical equipment and high frequency surgical accessories

201.1 Scope, object and related standards

Clause 1 of the general standard¹⁾ applies, except as follows:

201.1.1 * Scope

Replacement:

This International Standard applies to the BASIC SAFETY and ESSENTIAL PERFORMANCE of HF SURGICAL EQUIPMENT as defined in 201.3.222.

HF SURGICAL EQUIPMENT having a RATED OUTPUT POWER not exceeding 50 W (for example for micro-COAGULATION, or for use in dentistry or ophthalmology) is exempt from certain of the requirements of this particular standard. These exemptions are indicated in the relevant requirements.

201.1.2 Object

Replacement:

The object of this particular standard is to establish particular BASIC SAFETY and ESSENTIAL PERFORMANCE requirements for HF SURGICAL EQUIPMENT as defined in 201.3.222.

201.1.3 Collateral standards

Addition:

This particular standard refers to those applicable collateral standards that are listed in Clause 2 of the general standard and Clause 2 of this particular standard.

IEC 60601-1-2 and IEC 60601-1-8 apply as modified in Clauses 202 and 208 respectively. IEC 60601-1-3, IEC 60601-1-10 and IEC 60601-1-11²⁾ do not apply. All other published collateral standards in the IEC 60601-1 series apply as published.

201.1.4 Particular standards

Replacement:

In the IEC 60601 series, particular standards may modify, replace or delete requirements contained in the general standard and collateral standards as appropriate for the particular

¹⁾ The general standard is IEC 60601-1:2005, *Medical electrical equipment – Part 1: General requirements for basic safety and essential performance*.

²⁾ IEC 60601-1-11, *Medical electrical equipment – Part 1-11: General requirements for basic safety and essential performance – Collateral Standard: Requirements for medical electrical equipment and medical electrical systems used in the home healthcare environment* (in preparation).

ME EQUIPMENT under consideration, and may add other BASIC SAFETY and ESSENTIAL PERFORMANCE requirements.

A requirement of a particular standard takes priority over the general standard.

For brevity, IEC 60601-1 is referred to in this particular standard as the general standard. Collateral standards are referred to by their document number.

The numbering of clauses and subclauses of this particular standard corresponds to that of the general standard with the prefix "201" (e.g. 201.1 in this standard addresses the content of Clause 1 of the general standard) or applicable collateral standard with the prefix "20x" where x is the final digit(s) of the collateral standard document number (e.g. 202.4 in this particular standard addresses the content of Clause 4 of the 60601-1-2 collateral standard, 203.4 in this particular standard addresses the content of Clause 4 of the 60601-1-3 collateral standard, etc.). The changes to the text of the general standard are specified by the use of the following words:

"Replacement" means that the clause or subclause of the general standard or applicable collateral standard is replaced completely by the text of this particular standard.

"Addition" means that the text of this particular standard is additional to the requirements of the general standard or applicable collateral standard.

"Amendment" means that the clause or subclause of the general standard or applicable collateral standard is amended as indicated by the text of this particular standard.

Subclauses, figures or tables which are additional to those of the general standard are numbered starting from 201.101. However, due to the fact that definitions in the general standard are numbered 3.1 through 3.139, additional definitions in this standard are numbered beginning from 201.3.201. Additional annexes are lettered AA, BB, etc., and additional items aa), bb), etc.

Subclauses or figures which are additional to those of a collateral standard are numbered starting from 20x, where "x" is the number of the collateral standard, e.g. 202 for IEC 60601-1-2, 203 for IEC 60601-1-3, etc.

The term "this standard" is used to make reference to the general standard, any applicable collateral standards and this particular standard taken together.

Where there is no corresponding clause or subclause in this particular standard, the clause or subclause of the general standard or applicable collateral standard, although possibly not relevant, applies without modification; where it is intended that any part of the general standard or applicable collateral standard, although possibly relevant, is not to be applied, a statement to that effect is given in this particular standard.

201.2 Normative references

NOTE Informative references are listed in the bibliography beginning on page 79.

Clause 2 of the general standard applies, except as follows:

Replacement:

IEC 60601-1-2:2007, *Medical electrical equipment – Part 1-2: General requirements for basic safety and essential performance – Collateral standard: Electromagnetic compatibility – Requirements and tests*

Customer: Bradley Hyde- No. of User(s): 1 - Company: Latham and Watkins
Order No.: WS-2023-008299 - IMPORTANT: This file is copyright of IEC, Geneva, Switzerland. All rights reserved.
This file is subject to a licence agreement. Enquiries to Email: sales@iec.ch - Tel.: +41 22 919 02 11

IEC 60601-1-8:2006, *Medical electrical equipment – Part 1-8: General requirements for basic safety and essential performance – Collateral standard: General requirements, tests and guidance for alarm systems in medical electrical equipment and medical electrical systems*

Addition:

IEC 61000-4-3:2006, *Electromagnetic compatibility (EMC) – Part 4-3: Testing and measurement techniques – Radiated, radio-frequency electromagnetic field immunity test*

IEC 61000-4-6:2003, *Electromagnetic compatibility (EMC) – Part 4-6: Testing and measurement techniques – Immunity to conducted disturbances, induced by radio-frequency fields*

CISPR 11:2003, *Industrial, scientific and medical (ISM) radio-frequency equipment – Electromagnetic disturbance characteristics – Limits and methods of measurement*

201.3 Terms and definitions

For the purposes of this document, the terms and definitions given in IEC 60601-1:2005, apply, except as follows:

Replace NOTE 1 with the following:

NOTE 1 Where the terms "voltage" and "current" are used in this document, they mean the r.m.s. values of an alternating, direct or composite voltage or current averaged over 1 s unless stated otherwise.

Addition:

201.3.201

ACTIVE ACCESSORY

HF SURGICAL ACCESSORY intended for manipulation by the OPERATOR to produce surgical effects at the intended site on the PATIENT, generally comprising an ACTIVE HANDLE, the cord of an ACTIVE ACCESSORY, ACTIVE CONNECTOR and ACTIVE ELECTRODE

201.3.202

ACTIVE CONNECTOR

part of an ACTIVE ACCESSORY intended for connection to an ACTIVE OUTPUT TERMINAL, which may include additional terminals for connection of a FINGERSWITCH to a SWITCH SENSOR

201.3.203

ACTIVE ELECTRODE

part of an ACTIVE ACCESSORY extending from the ACTIVE HANDLE to the surgical site

201.3.204

ACTIVE ELECTRODE INSULATION

electrical insulation material affixed to part of an ACTIVE ELECTRODE intended to prevent unintended injury to the OPERATOR or adjacent PATIENT tissue

201.3.205

ACTIVE HANDLE

part of an ACTIVE ACCESSORY intended to be held by the OPERATOR

201.3.206

ACTIVE OUTPUT TERMINAL

part of HF SURGICAL EQUIPMENT or ASSOCIATED EQUIPMENT intended for connection to an ACTIVE ACCESSORY and for delivery of HF current thereto

Customer: Bradley Hyde- No. of User(s): 1 - Company: Latham and Watkins

Order No.: WS-2023-008299 - IMPORTANT: This file is copyright of IEC, Geneva, Switzerland. All rights reserved.

This file is subject to a licence agreement. Enquiries to Email: sales@iec.ch - Tel.: +41 22 919 02 11

201.3.207***ASSOCIATED EQUIPMENT**

equipment other than HF SURGICAL EQUIPMENT that may be electrically connected to the PATIENT circuit and not intended for independent use

201.3.208***BIPOLAR**

method of applying HF output current to a PATIENT via multiple-pole ACTIVE ELECTRODES

201.3.209**BIPOLAR ELECTRODE**

assembly of two or more ACTIVE ELECTRODES on the same support, so constructed that, when energized, the HF current flows mainly amongst these electrodes

201.3.210**COAGULATION**

use of HF current to elevate the temperature of tissue, e.g. to reduce or terminate undesired bleeding

NOTE COAGULATION may take the form of contact or non-contact COAGULATION.

201.3.211**CONTACT QUALITY MONITOR****CQM**

circuit in HF SURGICAL EQUIPMENT or ASSOCIATED EQUIPMENT intended for connection to a MONITORING NE providing an alarm in the event that NEUTRAL ELECTRODE (NE) contact with the PATIENT becomes insufficient

NOTE A CONTACT QUALITY MONITOR is functional only when used with a MONITORING NE.

201.3.212**CONTINUITY MONITOR**

circuit in HF SURGICAL EQUIPMENT or ASSOCIATED EQUIPMENT intended for connection to an NE, except MONITORING NE, providing an alarm in the event of electrical discontinuity in the NE cable or its connections

201.3.213***CREST FACTOR**

dimensionless value equal to the peak output voltage divided by the r.m.s. voltage as measured at the output of HF SURGICAL EQUIPMENT in an open circuit condition

NOTE Specific information on the correct way to make the measurements needed to calculate this value may be found in Annex AA.

201.3.214***CUTTING**

resection or dissection of body tissue caused by the passage of HIGH FREQUENCY current of high current density at the ACTIVE ELECTRODE(S)

201.3.215***EARTH REFERENCED PATIENT CIRCUIT**

PATIENT circuit which includes components, such as capacitors, installed to provide a low-impedance path to earth for HF currents

201.3.216**FINGERSWITCH**

device generally included with an ACTIVE ACCESSORY which, when manipulated by the OPERATOR, enables HF output to be produced and, when released disables HF output

Customer: Bradley Hyde- No. of User(s): 1 - Company: Latham and Watkins

Order No.: WS-2023-008299 - IMPORTANT: This file is copyright of IEC, Geneva, Switzerland. All rights reserved.

This file is subject to a licence agreement. Enquiries to Email: sales@iec.ch - Tel.: +41 22 919 02 11

SERENDIA-00084161

Ex. 2077-013

CX-2848.13

NOTE Requirements for similar switches intended to perform functions other than activation of HF output are under consideration.

201.3.217

*FULGURATION

form of COAGULATION using long (0,5 mm or more) electrical sparks to heat tissue surfaces superficially, with no intentional mechanical contact between the ACTIVE ELECTRODE and the tissue

201.3.218

*HIGH FREQUENCY

HF

frequencies generally greater than 200 kHz

201.3.219

HF ISOLATED PATIENT CIRCUIT

PATIENT circuit where there are no components installed to provide a low-impedance path to earth for HF currents

201.3.220

HF PATIENT CIRCUIT

any electrical circuit which contains one or more PATIENT CONNECTIONS

201.3.221

HF SURGICAL ACCESSORY

ACCESSORY intended to conduct, supplement or monitor HF energy applied to the PATIENT from HF SURGICAL EQUIPMENT

NOTE HF SURGICAL ACCESSORIES include HF surgical application electrodes, including cords and connectors for attachment to HF SURGICAL EQUIPMENT, as well as other ASSOCIATED EQUIPMENT intended for connection to the HF surgical PATIENT circuit.

201.3.222

HF SURGICAL EQUIPMENT

MEDICAL ELECTRICAL EQUIPMENT including associated ACCESSORIES intended for the performance of surgical operations, such as the CUTTING or COAGULATION of biological tissue by means of HIGH FREQUENCY currents

NOTE HF SURGICAL EQUIPMENT is also known as surgical diathermy or electrosurgical equipment.

201.3.223

*HF SURGICAL MODE

any of a number of OPERATOR selectable HF output characteristics intended to provide a specific indicated surgical effect at a connected ACTIVE ACCESSORY, such as CUTTING, COAGULATION and the like

NOTE Each available HF SURGICAL MODE may be provided with an OPERATOR-adjustable output control to set the desired intensity or speed of the surgical effect.

201.3.224

*MAXIMUM OUTPUT VOLTAGE

for each available HF SURGICAL MODE, the magnitude of the maximum possible peak HF output voltage appearing between PATIENT circuit connections

201.3.225

*MONITORING NE

NE intended for use with a CONTACT QUALITY MONITOR

Customer: Bradley Hyde- No. of User(s): 1 - Company: Latham and Watkins

Order No.: WS-2023-008299 - IMPORTANT: This file is copyright of IEC, Geneva, Switzerland. All rights reserved.

This file is subject to a licence agreement. Enquiries to Email: sales@iec.ch - Tel.: +41 22 919 02 11

201.3.226***MONOPOLAR**

method of applying HF output current to a PATIENT via an ACTIVE ELECTRODE and returning via a separately-connected NEUTRAL ELECTRODE or via the PATIENT'S body capacitance to earth

201.3.227**NEUTRAL ELECTRODE****NE**

electrode of a relatively large area for connection to the body of the PATIENT, intended to provide a return path for the HIGH FREQUENCY current with such a low current density in the body tissue that physical effects such as unwanted burns are avoided

NOTE The NEUTRAL ELECTRODE is also known as plate, plate electrode, passive, return or dispersive electrode.

201.3.228**RATED ACCESSORY VOLTAGE**

maximum peak HF output voltage which may be applied to a MONOPOLAR HF SURGICAL ACCESSORY with respect to an NE connected to the PATIENT. For a BIPOLAR HF SURGICAL ACCESSORY, the maximum peak HF output voltage which may be applied to pairs of opposite polarity

201.3.229**RATED LOAD**

value of non-reactive load resistance which, when connected, results in the maximum HF output power from each HF SURGICAL MODE of the HF SURGICAL EQUIPMENT

201.3.230**RATED OUTPUT POWER**

for each HF SURGICAL MODE set at its maximum output setting, the power in watts produced when all ACTIVE OUTPUT TERMINALS which can be activated simultaneously are connected to their respective RATED LOADS

201.3.231**SWITCH SENSOR**

part of HF SURGICAL EQUIPMENT or ASSOCIATED EQUIPMENT which controls activation of HF output in response to operation of a connected FINGERSWITCH or footswitch

201.4 General requirements

Clause 4 of the general standard applies, except as follows:

201.4.2 * RISK MANAGEMENT PROCESS for ME EQUIPMENT or ME SYSTEMS

Addition:

MANUFACTURERS shall include, within their RISK ANALYSIS, the potential for their HF SURGICAL EQUIPMENT and/or HF SURGICAL ACCESSORIES to be used in MONOPOLAR high current situations and the impact this would have on the heating under the NEUTRAL ELECTRODE (for example, see 201.7.9.2.2.101 f)).

201.4.3 * ESSENTIAL PERFORMANCE

Addition:

NOTE 101 Please refer to Annex AA.

201.4.7 SINGLE FAULT CONDITION for ME EQUIPMENT

Addition:

Customer: Bradley Hyde- No. of User(s): 1 - Company: Latham and Watkins

Order No.: WS-2023-008299 - IMPORTANT: This file is copyright of IEC, Geneva, Switzerland. All rights reserved.

This file is subject to a licence agreement. Enquiries to Email: sales@iec.ch - Tel.: +41 22 919 02 11

Additional SINGLE FAULT CONDITIONS:

- aa) failure in the CONTINUITY MONITOR or CONTACT QUALITY MONITOR which might cause a unacceptable RISK (see 201.8.4.101);
- bb) a defect in the output switching circuit resulting in an excessive low-frequency PATIENT LEAKAGE CURRENT (see 201.8.10.4.101.1);
- cc) any defect which results in the unwanted energization of the PATIENT circuit (see 201.12.4.2.101);
- dd) any defect which results in a significant increase in output power relative to the output setting (see 201.12.4.4.101).

201.4.11 Power input

Amendment:

The HF SURGICAL EQUIPMENT shall be operated in the output mode and using the load which creates the greatest steady state current.

201.5 General requirements for testing of ME EQUIPMENT

Clause 5 of the general standard applies, except as follows:

201.5.4 * Other conditions

Addition:

- aa) Particular care shall be taken to ensure accuracy and safety during measurement of HF output. See Annex AA for guidance.

201.6 Classification of ME EQUIPMENT and ME SYSTEMS

Clause 6 of the general standard applies.

201.7 ME EQUIPMENT identification, marking and documents

Clause 7 of the general standard applies, except as follows:

201.7.2.8.2 Other power sources

Amendment:

Subclause 7.2.8.2 of the general standard does not apply to ACTIVE OUTPUT TERMINALS.

201.7.2.10 APPLIED PARTS

Addition:

The relevant symbols required for marking DEFIBRILLATION-PROOF APPLIED PARTS shall be attached to the front panel, but are not required on the APPLIED PARTS.

Connections on the HF SURGICAL EQUIPMENT and ASSOCIATED EQUIPMENT for the connection of NE leads shall be marked with the following symbols:



IEC 1192/06

Figure 201.101 – Symbol used with an EARTH REFERENCED PATIENT CIRCUIT

Customer: Bradley Hyde- No. of User(s): 1 - Company: Latham and Watkins

Order No.: WS-2023-008299 - IMPORTANT: This file is copyright of IEC, Geneva, Switzerland. All rights reserved.

This file is subject to a licence agreement. Enquiries to Email: sales@iec.ch - Tel.: +41 22 919 02 11



IEC 1193/06

Figure 201.102 – Symbol used with a HF ISOLATED PATIENT CIRCUIT**201.7.4.2 * Control devices***Addition:*

The output control shall have a scale and/or associated indicator showing the relative units of HIGH FREQUENCY output. The indication shall not be marked in watts unless the indicated power is delivered with an accuracy of $\pm 20\%$ over the total load resistance range specified in 201.7.9.3.1.

The numeral "0" shall not be used unless no HF power in excess of 10 mW is delivered from an ACTIVE ELECTRODE or BIPOLAR ELECTRODE in this position.

NOTE The compliance test is the application of subclause 201.12.1.102.

201.7.8.1 * Colours of indicator lights*Replace Table 2 in the general standard with the following:***Table 201.101 – Colours of indicator lights and their meaning for HF SURGICAL EQUIPMENT**

Colour	Meaning
Red	Warning – immediate response by the OPERATOR is required, for example, a fault in the PATIENT circuit
Yellow	CUTTING mode
Blue	COAGULATION mode
Green	Ready for use
Any other colour	Meaning other than that of red, yellow, blue or green

201.7.8.2 * Colours of controls*Addition:*

Where operating controls, output terminals, indicator lights, pedals (see 201.12.2) and pushbuttons of FINGERSWITCHES (see 201.12.2) are associated with a particular HF SURGICAL MODE, they shall be identified by a consistent, unique colour not in conflict with Table 201.101.

Compliance is checked by inspection.

201.7.9.2.2 Warning and safety notices*Additional subclause:***201.7.9.2.2.101 * Additional information in instructions for use**

a) Notes on the application of HF SURGICAL EQUIPMENT. These notes shall draw the attention of the OPERATOR to certain precautions which are necessary in order to reduce the risk of accidental burns. In particular, advice, when appropriate, shall be given on the following:

- 1) * The entire area of the NEUTRAL ELECTRODE should be reliably attached to a suitably prepared and appropriate area of the PATIENT'S body as defined by the MANUFACTURER.

Customer: Bradley Hyde- No. of User(s): 1 - Company: Latham and Watkins

Order No.: WS-2023-008299 - IMPORTANT: This file is copyright of IEC, Geneva, Switzerland. All rights reserved.

This file is subject to a licence agreement. Enquiries to Email: sales@iec.ch - Tel.: +41 22 919 02 11

- 2) * The PATIENT should not come into contact with metal parts which are earthed or which have an appreciable capacitance to earth (for example operating table supports, etc.).
 - 3) * Skin-to-skin contact (for example between the arms and body of the PATIENT) should be avoided, for example by insertion of dry gauze.
 - 4) * When HF SURGICAL EQUIPMENT and physiological monitoring equipment are used simultaneously on the same PATIENT, any monitoring electrodes should be placed as far as possible from the surgical electrodes. Needle monitoring electrodes are not recommended.
In all cases, monitoring systems incorporating HIGH FREQUENCY current limiting devices are recommended.
 - 5) * The PATIENT leads should be positioned in such a way that contact with the PATIENT or other leads is avoided.
Temporarily unused ACTIVE ELECTRODES should be stored in a location that is isolated from the PATIENT.
 - 6) * For surgical procedures where the HF current could flow through parts of the body having a relatively small cross sectional area, the use of BIPOLAR techniques may be desirable in order to avoid unwanted tissue damage.
 - 7) The output power selected should be as low as possible for the intended purpose. Certain devices or ACCESSORIES may present an unacceptable RISK at low power settings. For example, with argon beam COAGULATION, the risk of gas embolism rises if there is insufficient HF power to produce a rapid, impermeable eschar on the target tissue.
 - 8) * Apparent low output or failure of the HF SURGICAL EQUIPMENT to function correctly at the normal operating settings may indicate faulty application of the NEUTRAL ELECTRODE or poor contact in its connections. In this case, the application of the NEUTRAL ELECTRODE and its connections should be checked before selecting a higher output power.
 - 9) The use of flammable anaesthetics or oxidizing gases such as nitrous oxide (N₂O) and oxygen should be avoided if a surgical procedure is carried out in the region of the thorax or the head, unless these agents are sucked away.
Non-flammable agents should be used for cleaning and disinfection wherever possible.
Flammable agents used for cleaning or disinfecting, or as solvents of adhesives, should be allowed to evaporate before the application of HF surgery. There is a risk of pooling of flammable solutions under the PATIENT or in body depressions such as the umbilicus, and in body cavities such as the vagina. Any fluid pooled in these areas should be mopped up before HF SURGICAL EQUIPMENT is used. Attention should be called to the danger of ignition of endogenous gases. Some materials, for example cotton and gauze, when saturated with oxygen may be ignited by sparks produced in NORMAL USE of the HF SURGICAL EQUIPMENT.
 - 10) For PATIENTS with cardiac pacemakers or other active implants, a possible hazard exists because interference with the action of the pacemaker may occur, or the pacemaker may be damaged. In case of doubt, approved qualified advice should be obtained.
 - 11) For HF SURGICAL EQUIPMENT with an operating mode as described in 201.12.2 c), a warning is required to the effect that the output from either ACTIVE ELECTRODE may change during use.
- b) A warning that interference produced by the operation of HF SURGICAL EQUIPMENT may adversely influence the operation of other electronic equipment.
- c) * For HF SURGICAL EQUIPMENT, the MAXIMUM OUTPUT VOLTAGE for each HF SURGICAL MODE and instruction regarding the RATED ACCESSORY VOLTAGE as follows:

- 1) For situations where the MAXIMUM OUTPUT VOLTAGE ($U_{\max}$) is less than or equal to 1 600 V, provide instruction that ASSOCIATED EQUIPMENT and ACTIVE ACCESSORIES should be selected that have a RATED ACCESSORY VOLTAGE equal to or greater than the MAXIMUM OUTPUT VOLTAGE.
- 2) For situations where the MAXIMUM OUTPUT VOLTAGE ($U_{\max}$) is greater than 1 600 V, calculate the variable y using the formula:

$$y = \frac{U_{\max} - 400 [\text{Volts}]}{600 [\text{Volts}]}$$

Take the smaller of variable y or the number 6. If the result is less than or equal to the CREST FACTOR for that HF SURGICAL MODE, then provide instruction that ASSOCIATED EQUIPMENT and ACTIVE ACCESSORIES should be selected that have a RATED ACCESSORY VOLTAGE equal to or greater than the MAXIMUM OUTPUT VOLTAGE.

- 3) For situations where the MAXIMUM OUTPUT VOLTAGE ($U_{\max}$) is greater than 1 600 V, and the CREST FACTOR is less than the variable y calculated above, a warning shall be provided that any ASSOCIATED EQUIPMENT and ACTIVE ACCESSORIES used with such mode or setting must be rated to withstand the combination of actual voltage and CREST FACTOR.

Where the MAXIMUM OUTPUT VOLTAGE varies with the output setting, that information shall be presented diagrammatically as a function of output setting.

- d) A warning that failure of the HF SURGICAL EQUIPMENT could result in an unintended increase of output power.
- e) * A statement of compatibility with specific MONITORING NE.
A warning that, unless a compatible MONITORING NE is used with a CONTACT QUALITY MONITOR, loss of safe contact between the NE and the PATIENT will not result in an auditory alarm.
NOTE 1 This requirement does not apply for HF SURGICAL EQUIPMENT only incorporating BIPOLAR output.
NOTE 2 This requirement does not apply for HF SURGICAL EQUIPMENT intended for use without a NEUTRAL ELECTRODE. (See 201.15.101).
- f) Where the temperature under the NEUTRAL ELECTRODE, during intended or foreseen use, may exceed the limits listed in subclause 11.1.2.2 of the general standard or 201.15.101.5 of this standard, instructions, warnings and cautions for proper use of the NEUTRAL ELECTRODE shall be provided.
- g) * A warning addressing the risks resulting from neuromuscular stimulation which can occur especially with modes which produce electrical arcs between the ACTIVE ELECTRODE and tissue.
- h) * For HF SURGICAL EQUIPMENT that can be energized without continuous activation of a SWITCH SENSOR as per subclause 201.8.10.4.101.2, warnings and/or cautions regarding the risks

201.7.9.2.14 * Accessories, supplementary equipment, used material

Addition:

The instructions for use shall include:

- a) Information concerning the selection and use of HF SURGICAL ACCESSORIES in order to avoid incompatibility and unsafe operation (see also 201.15.4.1.101 and 201.15.4.1.102).
- b) Advice for the OPERATOR to avoid HF output settings where MAXIMUM OUTPUT VOLTAGE may exceed RATED ACCESSORY VOLTAGE.
- c) Advice concerning the compatibility between a MONITORING NE and a CONTACT QUALITY MONITOR.
- d) Advice for the OPERATOR regularly to inspect the ACCESSORIES. In particular, electrode cables and HF ENERGIZED ENDOTHERAPY DEVICES (see IEC 60601-2-18) should be checked (e.g. under magnification) for possible damage.

Customer: Bradley Hyde- No. of User(s): 1 - Company: Latham and Watkins

Order No.: WS-2023-008299 - IMPORTANT: This file is copyright of IEC, Geneva, Switzerland. All rights reserved.

This file is subject to a licence agreement. Enquiries to Email: sales@iec.ch - Tel.: +41 22 919 02 11

- e) * For ASSOCIATED EQUIPMENT and ACTIVE ACCESSORIES, including separately supplied parts thereof, the RATED ACCESSORY VOLTAGE .
- f) * On end use packaging for NEUTRAL ELECTRODES:
 - If marked for single use, an expiration date.
 - Information necessary to prevent burns at the NE site, e.g. limitation of output setting, PATIENT preparation and/or activation duration.
 - If intended for use only on small PATIENTS, a marking in kg indicating the maximum PATIENT weight for which it is intended to be used. See 201.15.101.5
- g) * On instructions for use for MONITORING NEUTRAL ELECTRODES:
 - A statement of compatibility with specific CONTACT QUALITY MONITOR(S).
- h) HF SURGICAL ACCESSORIES where the temperature under the NE, during intended or foreseen use, may result in the temperature exceeding the limits listed in subclause 11.1.2.2 of the general standard or subclause 201.15.101.5 of this standard shall be accompanied by instructions, warnings and cautions for the proper use of NEUTRAL ELECTRODES.
- i) On instructions for use for HF SURGICAL ACCESSORIES intended for use only with specific HF SURGICAL EQUIPMENT or HF waveforms or voltages, a detailed statement to that effect.

201.7.9.2.15 Environmental protection

Addition:

- Provide advice to the OPERATOR regarding the advisability of the use of smoke-plume extraction.

201.7.9.3 Technical description

201.7.9.3.1 * General

Addition:

- power output data – MONOPOLAR output (for all HF SURGICAL MODES available, any variable "blend" control being set to the maximum position) including:
 - diagrams showing the power output at full and half output control settings minimally over the range of load resistance 100 Ω to 2 000 Ω , but extended as necessary to include the RATED LOAD;
 - diagrams showing the power output versus the output control setting at a specified load resistance in the range as defined above;
- power output data – BIPOLAR output (for all HF SURGICAL MODES as defined above) including:
 - diagrams showing the power output at full and half output control settings minimally over the range of load resistance 10 Ω to 1 000 Ω , but extended as necessary to include the RATED LOAD;
 - diagrams showing the power output versus the output control setting at a specified load resistance in the range as defined above;
- voltage output data – MONOPOLAR and BIPOLAR output (for all HF SURGICAL MODES available). Maximum voltage data required by 201.7.9.2.2.101 c);
- where HF SURGICAL EQUIPMENT is specified for use without a NEUTRAL ELECTRODE, this shall be stated;
- where HF SURGICAL EQUIPMENT or ASSOCIATED EQUIPMENT is designed to have a single, fixed output setting, then reference to "half output control settings" shall be ignored.

Customer: Bradley Hyde- No. of User(s): 1 - Company: Latham and Watkins

Order No.: WS-2023-008299 - IMPORTANT: This file is copyright of IEC, Geneva, Switzerland. All rights reserved.

This file is subject to a licence agreement. Enquiries to Email: sales@iec.ch - Tel.: +41 22 919 02 11

SERENDIA-00084168

Ex. 2077-020

CX-2848.20

201.8 Protection against electrical HAZARDS from ME EQUIPMENT

Clause 8 of the general standard applies, except as follows:

201.8.3 Classification of APPLIED PARTS

Addition:

- aa) The APPLIED PARTS of HF SURGICAL EQUIPMENT shall be TYPE BF or TYPE CF APPLIED PARTS.

Additional subclause:

201.8.4.101 * NEUTRAL ELECTRODE monitoring circuit

HF SURGICAL EQUIPMENT having a RATED OUTPUT POWER of more than 50 W shall be provided with a CONTINUITY MONITOR and/or a CONTACT QUALITY MONITOR arranged so as to de-energize the output and to give an audible alarm when a failure of the NEUTRAL ELECTRODE circuit or its connections occurs. The audible alarm shall meet the sound level requirements of 201.12.4.2.101 and shall not be externally adjustable.

NOTE This audible alarm and visible indicator light are not intended to meet the definition of an ALARM SIGNAL in IEC 60601-1-8. See also Clause 208 of this standard.

The monitoring circuit shall be supplied from a power source isolated from the MAINS PART and from earth and having a voltage not exceeding 12 V. The limitation of monitoring current for a CONTACT QUALITY MONITOR is defined in 201.8.7.3.

An additional visible warning consisting of a red indicator light shall be provided (see 201.7.8.1).

Compliance of a CONTINUITY MONITOR is checked by operating the HF SURGICAL EQUIPMENT at maximum output control setting in each operating mode into the circuit shown in Figure 201.103. The switch is closed and opened five times and the HF output shall be disabled and the alarm shall sound at each opening of the switch.

Compliance of a CONTACT QUALITY MONITOR is checked by switching on the mains of the HF SURGICAL EQUIPMENT and setting its controls for MONOPOLAR operation, except that it shall not be activated. Then a compatible MONITORING NE, selected according to the advice per 201.7.9.2.2.101 e), is connected to the NE connections of the CONTACT QUALITY MONITOR. The NE is then placed, according to marked instructions for use, with full contact on a human subject or a suitable surrogate surface, and the CONTACT QUALITY MONITOR is set up according to instructions for use. The HF SURGICAL EQUIPMENT is then activated in a MONOPOLAR HF SURGICAL MODE. No alarm shall sound and HF output shall be present. With the HF SURGICAL EQUIPMENT now activated, the contact area between the NE and the human subject or a suitable surrogate surface is gradually reduced until a NE alarm occurs. The remaining contact area (alarm area), A_a shall be recorded for subsequent thermal rise testing per subclause 201.15.101.5, and no HF output shall be produced when activation is attempted.

This test shall be repeated along both axes using at least three samples of each compatible MONITORING NE.

Customer: Bradley Hyde- No. of User(s): 1 - Company: Latham and Watkins


Order No.: WS-2023-008299 - IMPORTANT: This file is copyright of IEC, Geneva, Switzerland. All rights reserved.

This file is subject to a licence agreement. Enquiries to Email: sales@iec.ch - Tel.: +41 22 919 02 11

SERENDIA-00084169

Ex. 2077-021

CX-2848.21

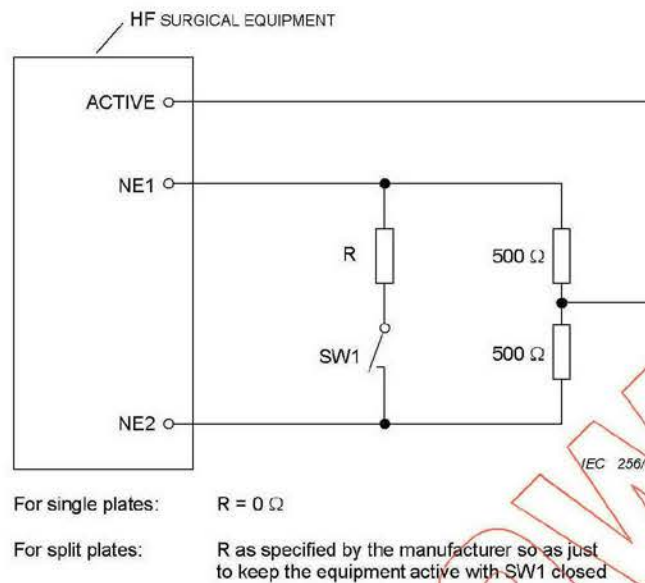

NOTE NEUTRAL ELECTRODES which are split into more than two parts should be tested accordingly.

Figure 201.103 – Circuit suitable for testing compliance to 201.8.4.101

201.8.4.102 * Neuromuscular stimulation

In order to minimize the possibility of neuromuscular stimulation, a capacitance shall be incorporated into the PATIENT circuit so that it is effectively in series with the ACTIVE ELECTRODE or one conductor of a BIPOLAR ELECTRODE. This capacitance shall not exceed 5 nF for MONOPOLAR PATIENT circuits and 50 nF for BIPOLAR PATIENT circuits. The d.c. resistance between ACTIVE and NEUTRAL ELECTRODE terminals, or between the terminals of a BIPOLAR output circuit, shall not be less than 2 MΩ.

Compliance is checked by inspection of the circuit arrangement and by measurement of the d.c. resistance between the output terminals.

201.8.5.1.2 * MEANS OF PATIENT PROTECTION (MOPP)

Amendment:

For HF SURGICAL EQUIPMENT and HF SURGICAL ACCESSORIES, the CREEPAGE DISTANCES and AIR CLEARANCES of insulation between the HF APPLIED PARTS and the ENCLOSURE including SIGNAL INPUT PARTS and SIGNAL OUTPUT PARTS, between the HF APPLIED PARTS and any intermediate circuit and between different HF PATIENT CIRCUITS shall be at least 3 mm/kV or 4 mm, whichever is the greater. The reference voltage shall be the maximum peak voltage. These separations need not be subjected to the dielectric strength test of 201.8.8.3.

NOTE: These creepage distances and air clearances are intended to represent two means of protection.

This requirement does not apply for components when the adequacy of ratings can be demonstrated, for example by component manufacturers' ratings or by the dielectric strength test of 201.8.8.3.

Customer: Bradley Hyde- No. of User(s): 1 - Company: Latham and Watkins
 Order No.: WS-2023-008299 - IMPORTANT: This file is copyright of IEC, Geneva, Switzerland. All rights reserved.
 This file is subject to a licence agreement. Enquiries to Email: sales@iec.ch - Tel.: +41 22 919 02 11

201.8.5.2.3 * PATIENT leads

Amendment:

This requirement shall not apply to the ACTIVE CONNECTORS or to any NE connectors except as detailed below.

For NEUTRAL ELECTRODE cables, the connector which is remote from the PATIENT shall be constructed so that the connections can not contact conductive live parts of FIXED mains socket outlets or MAINS CONNECTORS.

If able to be plugged into a FIXED mains socket-outlet or MAINS CONNECTOR, the said part shall be protected from making contact with parts at mains voltage by insulating means providing a CREEPAGE DISTANCE of at least 1,0 mm and a dielectric strength of 1 500 V.

Compliance is checked by inspection and by applying the dielectric strength test to the conductive connection of that part of the connector identified above.

201.8.5.5 DEFIBRILLATION-PROOF APPLIED PARTS

Amendment:

PATIENT circuits of HF SURGICAL EQUIPMENT shall be considered as APPLIED PARTS in the context of this subclause.

Compliance is checked by the common-mode test only, as described in subclause 8.5.5.1 and Figure 9 of the general standard using a test voltage of 2 kV instead of 5 kV.

After this test, HF SURGICAL EQUIPMENT shall be capable of meeting all the requirements and tests of this standard and of performing its intended function as described in the ACCOMPANYING DOCUMENTS.

201.8.6.1 * Applicability of requirements

Addition:

Generally, a PROTECTIVE EARTH CONDUCTOR shall not carry functional current. However, in HF SURGICAL EQUIPMENT having a RATED OUTPUT POWER not exceeding 50 W and intended for use without a NEUTRAL ELECTRODE, the PROTECTIVE EARTH CONDUCTOR of the mains cord may be used as a return path for the functional HIGH FREQUENCY current.

201.8.7.1 * General requirements

Item b)

Addition:

- with the HF not energized, but in such a way that the low-frequency LEAKAGE CURRENTS are not affected.

Amendment:

These investigations shall be carried out with the HF SURGICAL EQUIPMENT switched on but with PATIENT circuits not activated.

201.8.7.3 * Allowable values

Item b)

Customer: Bradley Hyde- No. of User(s): 1 - Company: Latham and Watkins

Order No.: WS-2023-008299 - IMPORTANT: This file is copyright of IEC, Geneva, Switzerland. All rights reserved.

This file is subject to a licence agreement. Enquiries to Email: sales@iec.ch - Tel.: +41 22 919 02 11

Addition:

PATIENT AUXILIARY CURRENTS associated with CONTACT QUALITY MONITORS shall not exceed the allowable values for TYPE BF APPLIED PARTS.

Item e)

Amendment:

The 10 mA limit for LEAKAGE CURRENT does not apply to HF LEAKAGE CURRENTS tested from ACTIVE and NEUTRAL ELECTRODES with PATIENT circuits activated (see 201.8.7.3.101).

Additional subclause:

201.8.7.3.101 Thermal effects of HF LEAKAGE CURRENTS

In order to prevent unintended thermal burns, HF LEAKAGE CURRENTS tested from ACTIVE and NEUTRAL ELECTRODES with HF PATIENT circuits activated shall, depending on their design, comply with the following requirements.

NOTE Independent requirements and conformance tests for HF SURGICAL EQUIPMENT and ACCESSORIES are under consideration.

***a) HIGH FREQUENCY LEAKAGE CURRENTS**


1) NEUTRAL ELECTRODE referenced to earth

The PATIENT circuit is isolated from earth but the NEUTRAL ELECTRODE is referenced to earth at HIGH FREQUENCIES by components (for example a capacitor) satisfying the requirements of a TYPE BF APPLIED PART. When tested as described below, the HF LEAKAGE CURRENT flowing from the NEUTRAL ELECTRODE through a non-inductive 200 Ω resistor to earth shall not exceed 150 mA.

Compliance is checked by the following tests.

Test 1 – The test is performed on each single output of the HF SURGICAL EQUIPMENT in turn with the electrode cables and electrodes as shown in Figure 201.104. The cables are spaced 0,5 m apart on an insulating surface 1 m above an earthed conductive plane.

The output is loaded with 200 Ω and the HF SURGICAL EQUIPMENT is operated at maximum output setting in each operating mode. The HF LEAKAGE CURRENT flowing from the NEUTRAL ELECTRODE through a non-inductive resistor of 200 Ω to earth is measured.

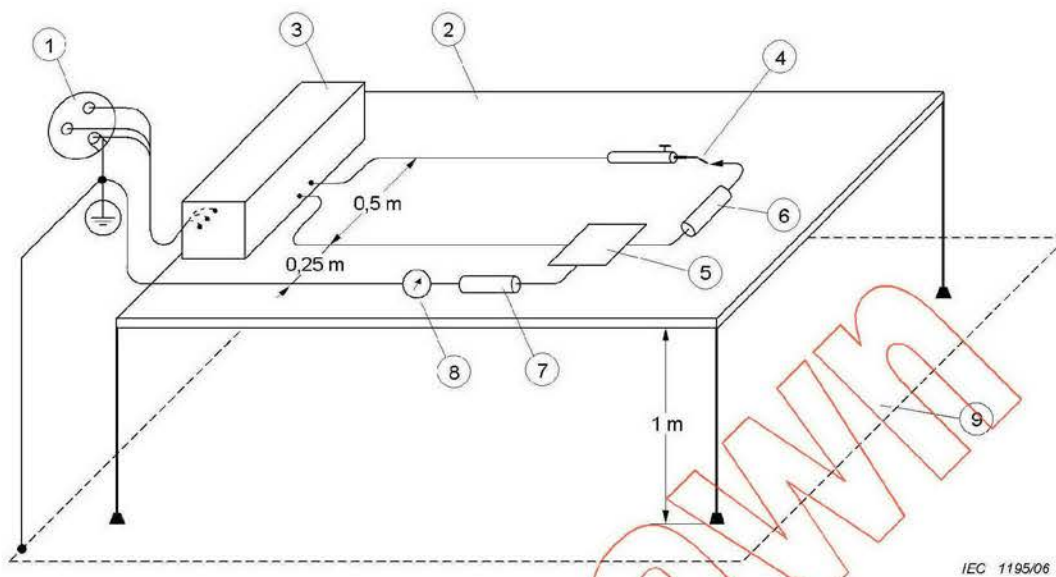


Key

- ① SUPPLY MAINS
- ② Table, made of insulating material
- ③ HF SURGICAL EQUIPMENT
- ④ ACTIVE ELECTRODE
- ⑤ NEUTRAL ELECTRODE, metallic or in contact with metal foil of the same size
- ⑥ Load resistance, 200 Ω
- ⑦ Measuring resistance, 200 Ω
- ⑧ HF current meter
- ⑨ Earthed conductive plane

Figure 201.104 – Measurement of HF LEAKAGE CURRENT with NEUTRAL ELECTRODE referenced to earth and load between electrodes

Test 2 – The HF SURGICAL EQUIPMENT is set up as for test 1, but the 200 Ω load resistor is connected between the ACTIVE ELECTRODE and the PROTECTIVE EARTH TERMINAL of the HF SURGICAL EQUIPMENT as shown in Figure 201.105. The HF LEAKAGE CURRENT flowing from the NEUTRAL ELECTRODE is measured.

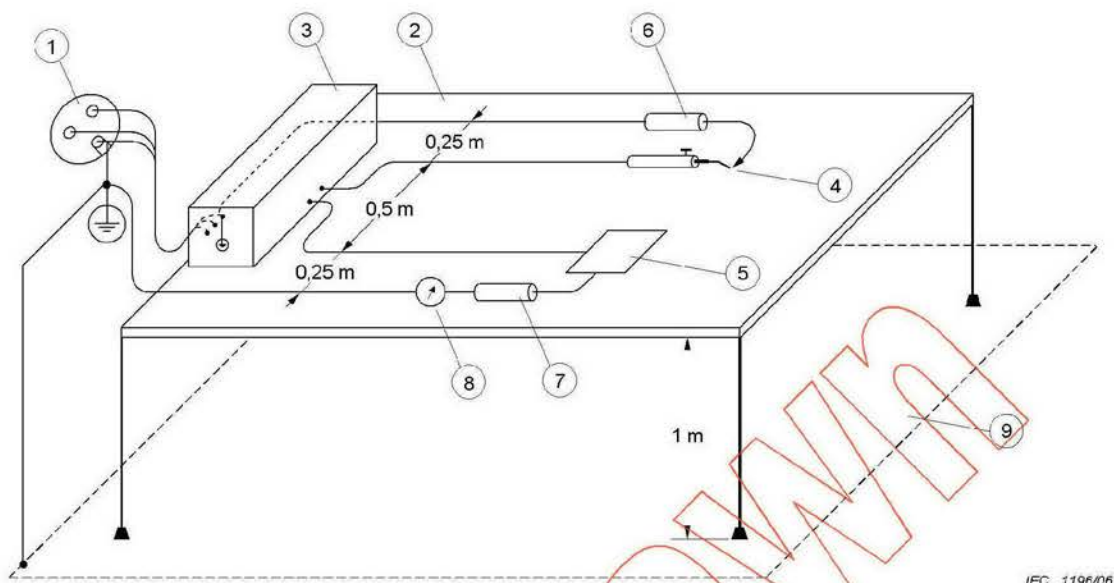

IEC 1196/06

Key

- ① SUPPLY MAINS
- ② Table, made of insulating material
- ③ HF SURGICAL EQUIPMENT
- ④ ACTIVE ELECTRODE
- ⑤ NEUTRAL ELECTRODE, metallic or in contact with metal foil of the same size
- ⑥ Load resistance, 200 Ω
- ⑦ Measuring resistance, 200 Ω
- ⑧ HF current meter
- ⑨ Earthed conductive plane

Figure 201.105 – Measurement of HF LEAKAGE CURRENT with NEUTRAL ELECTRODE referenced to earth and load from ACTIVE ELECTRODE to earth

2) NEUTRAL ELECTRODE isolated from earth at HIGH FREQUENCY


The PATIENT circuit is isolated from earth at both high and low frequencies, and the isolation shall be such that the HF LEAKAGE CURRENT flowing, in turn, from each electrode through a 200 Ω non-inductive resistor to earth does not exceed 150 mA when tested as described below.

Compliance is checked by the following test.

The HF SURGICAL EQUIPMENT is set up as shown in Figure 201.106, the output being unloaded and loaded at the RATED LOAD.

Any metal ENCLOSURES of CLASS II HF SURGICAL EQUIPMENT and INTERNALLY POWERED HF SURGICAL EQUIPMENT shall be connected to earth. HF SURGICAL EQUIPMENT having an insulating ENCLOSURE shall be positioned on earthed metal having an area at least equal to that of the base of the HF SURGICAL EQUIPMENT, during this test. The HF LEAKAGE CURRENT is measured from each electrode in turn while the HF SURGICAL EQUIPMENT is operated at maximum output setting in each HF SURGICAL MODE.

NOTE The above requirements do not apply for HF SURGICAL EQUIPMENT having a RATED OUTPUT POWER not exceeding 50 W and intended for use without a NEUTRAL ELECTRODE.

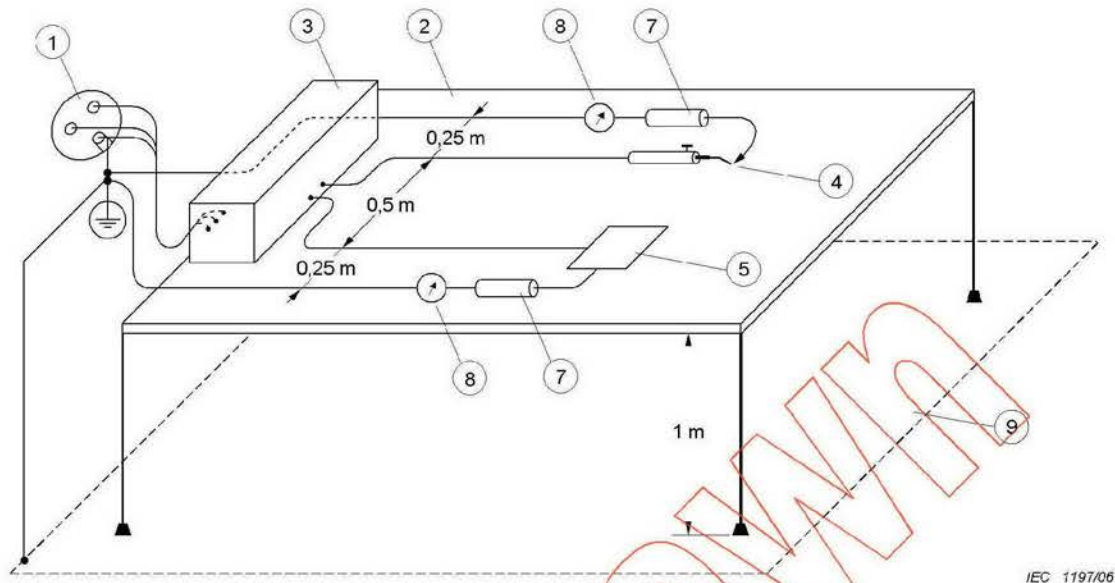

IEC 1197/06

Key

- ① SUPPLY MAINS
- ② Table, made of insulating material
- ③ HF SURGICAL EQUIPMENT
- ④ ACTIVE ELECTRODE
- ⑤ NEUTRAL ELECTRODE, metallic or in contact with metal foil of the same size
- ⑦ Measuring resistance, 200 Ω
- ⑧ HF current meter
- ⑨ Earthed conductive plane

Figure 201.106 – Measurement of HF LEAKAGE CURRENT with NEUTRAL ELECTRODE isolated from earth at HIGH FREQUENCY

***3) BIPOLAR application**

Any PATIENT circuit specifically designed for BIPOLAR application shall be isolated from earth and from other APPLIED PARTS at both high and low frequencies.

The HF LEAKAGE CURRENT flowing from either pole of the BIPOLAR output to earth and to the NEUTRAL ELECTRODE via a 200 Ω non-inductive resistor in each line shall not exceed the value which produces a power in a 200 Ω non-inductive resistor equal to 1 % of the maximum BIPOLAR RATED OUTPUT POWER, with all output controls set to maximum.

Compliance is checked by the following test.

The HF SURGICAL EQUIPMENT is set up as shown in Figure 201.107. The test is conducted using one side of the BIPOLAR output and using BIPOLAR and (if applicable) NEUTRAL ELECTRODE leads supplied or recommended by the manufacturer. The test is conducted with the output first being unloaded and then repeated with the output loaded at the RATED LOAD. The squared current value multiplied by 200 Ω shall not exceed the requirement above. The test is then repeated for the other side of the BIPOLAR output.

Any metal ENCLOSURES of CLASS II HF SURGICAL EQUIPMENT and INTERNALLY POWERED HF SURGICAL EQUIPMENT shall be connected to earth. HF SURGICAL EQUIPMENT having an insulating ENCLOSURE shall be positioned on earthed metal having an area at least equal to that of the base of the HF SURGICAL EQUIPMENT.

Customer: Bradley Hyde- No. of User(s): 1 - Company: Latham and Watkins

Order No.: WS-2023-008299 - IMPORTANT: This file is copyright of IEC, Geneva, Switzerland. All rights reserved.

This file is subject to a licence agreement. Enquiries to Email: sales@iec.ch - Tel.: +41 22 919 02 11


SERENDIA-00084175

Ex. 2077-027

CX-2848.27

During all measurements of HF LEAKAGE CURRENTS, the POWER SUPPLY CORD of the HF SURGICAL EQUIPMENT shall be folded up to form a bundle having a length not exceeding 40 cm.

NOTE The above requirements 1), 2) and 3) apply to HF SURGICAL EQUIPMENT with both TYPE BF and TYPE CF APPLIED PARTS. Requirements for HF ENCLOSURE LEAKAGE CURRENTS are under consideration.

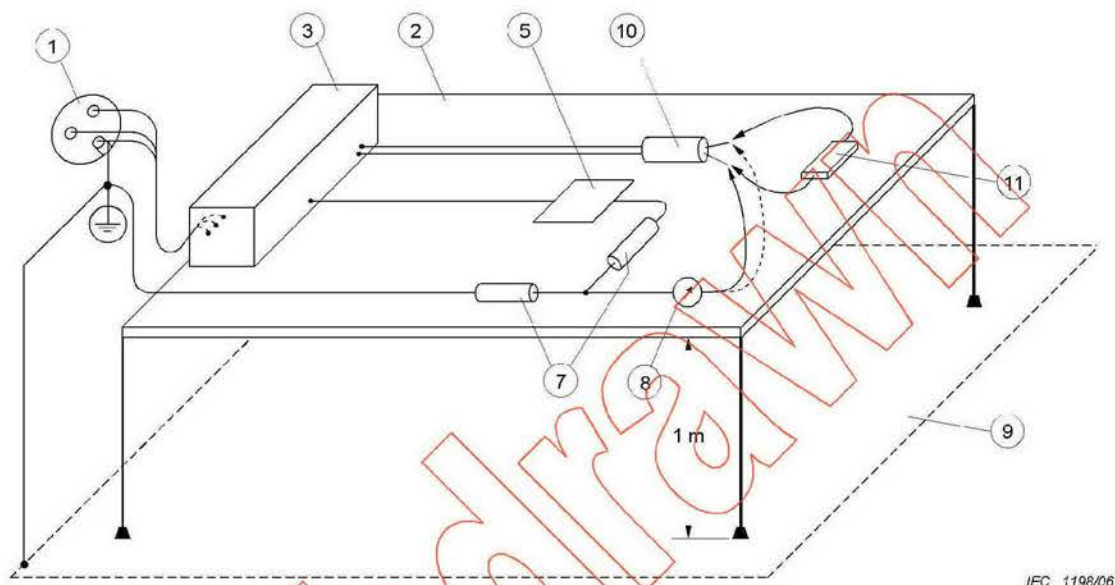

Key

- ① SUPPLY MAINS
- ② Table, made of insulating material
- ③ HF SURGICAL EQUIPMENT
- ⑤ NEUTRAL ELECTRODE, metallic or in contact with metal foil of the same size
- ⑦ Measuring resistance, 200 Ω
- ⑧ HF current meter
- ⑨ Earthed conductive plane
- ⑩ Activated BIPOLAR ELECTRODE
- ⑪ Load resistance as required with HF power measuring device

Figure 201.107 – Measurement of HF LEAKAGE CURRENT from a BIPOLAR ELECTRODE

*b) HIGH FREQUENCY LEAKAGE CURRENTS measured directly at the HF SURGICAL EQUIPMENT terminals

The preceding item a) may alternatively be fulfilled with a limit of 100 mA for 1) and 2) and with unchanged limits corresponding to 1 % of the BIPOLAR RATED OUTPUT POWER into 200 Ω and not to exceed 100 mA for 3) when the HF LEAKAGE CURRENT is measured directly at the HF SURGICAL EQUIPMENT terminals.

Compliance is checked by measurement similar to the tests described in 201.8.7.3.101 a), but without the electrode cables, and using leads as short as practicable for connecting the load resistor, the measuring resistor and the current measuring instrument to the HF SURGICAL EQUIPMENT terminals.

c) Cross-coupling between different HF PATIENT CIRCUITS

Customer: Bradley Hyde- No. of User(s): 1 - Company: Latham and Watkins

Order No.: WS-2023-008299 - IMPORTANT: This file is copyright of IEC, Geneva, Switzerland. All rights reserved.

This file is subject to a licence agreement. Enquiries to Email: sales@iec.ch - Tel.: +41 22 919 02 11

When any other PATIENT circuit is activated at its highest output settings and at all available operation modes, then:

- 1) A non-activated MONOPOLAR PATIENT circuit shall produce no more than 150 mA HIGH FREQUENCY current into a 200 Ω load to earth and, in turn, to the NEUTRAL ELECTRODE.
- 2) A non-activated BIPOLAR PATIENT circuit shall produce no more than 50 mA into a 200 Ω load connected across the two terminals or – with short circuited terminals – into a 200 Ω load to earth and into a 200 Ω load to the NEUTRAL ELECTRODE (both currents added, see Figure 201.107).

Compliance is checked by measurements using the test arrangements specified in subclause 201.8.7.3.101 b) and the HF SURGICAL EQUIPMENT is set up as shown in Figure 201.106 (for MONOPOLAR) or Figure 201.107 (for BIPOLAR PATIENT circuits).

201.8.8.3 Dielectric strength

Amendment:

The requirements and tests for HF SURGICAL ACCESORIES are given in 201.8.8.3.101 and 201.15.101.4.

Additional test conditions:

- aa) If, during dielectric strength testing of solid insulation forming MEANS OF PATIENT PROTECTION, a breakdown or flashover occurs through the atmosphere at the AIR CLEARANCE specified in 8.9 of the general standard and 201.8.5.1.2 of this standard, an insulating barrier may be placed to prevent this breakdown so that the protective insulation can be tested.
- bb) If, during dielectric strength testing of solid insulation forming MEANS OF PATIENT PROTECTION, a breakdown or flashover occurs at the CREEPAGE DISTANCE specified in 8.9 of the general standard and 201.8.5.1.2 of this standard, then the test shall be carried out on the components which provide MEANS OF PATIENT PROTECTION, such as transformers, relays, optocouplers or CREEPAGE DISTANCES on printed circuit boards.

Additional subclauses:

201.8.8.3.101 * ACTIVE ACCESSORY insulation

ACTIVE ACCESSORIES and cords of ACTIVE ACCESSORIES shall be sufficiently insulated to mitigate unintended thermal burn risk to the PATIENT and OPERATOR under conditions of NORMAL USE.

Compliance is checked as follows:

Test samples, other than those marked for single use, shall have undergone the sterilization methods and number of cycles as outlined in the instructions for use. See 7.9.2.12 in the general standard.

The insulated parts of all ACTIVE ACCESSORIES, other than ACTIVE HANDLES and ACTIVE CONNECTORS, shall be preconditioned by immersion in 0,9 % saline for 12 h. Operative conductors which may have been exposed in preparation for testing, as well as the insulation of the cords of ACTIVE ACCESSORIES within 100 mm of the ends, shall be protected from contact with saline. Upon completion of this preconditioning, excess saline shall be removed from surfaces and cavities by shaking and/or wiping with a dry cloth.

Immediately following saline preconditioning, applicable electrical testing shall be conducted in the following order:

- HF leakage (201.8.8.3.102);

Customer: Bradley Hyde- No. of User(s): 1 - Company: Latham and Watkins

Order No.: WS-2023-008299 - IMPORTANT: This file is copyright of IEC, Geneva, Switzerland. All rights reserved.

This file is subject to a licence agreement. Enquiries to Email: sales@iec.ch - Tel.: +41 22 919 02 11

- HF dielectric strength (201.8.8.3.103);
- mains frequency dielectric strength (201.8.8.3.104).

201.8.8.3.102 * ACTIVE ACCESSORY HF leakage

a) Measured HF LEAKAGE CURRENT

The insulation applied to cords for ACTIVE ACCESSORIES intended for MONOPOLAR application shall limit HF LEAKAGE CURRENT I_{leakage} passing through the external surface of the insulation to less than

$$I_{\text{leakage}} [\text{mA}] = 9,0 \times 10^{-6} \times d \times L \times f_{\text{test}} \times U_{\text{peak}} \text{ where}$$

d = smallest outer dimension of the insulation in mm,

f_{test} = HF test voltage frequency in kHz,

L = length of sample insulation through which HF LEAKAGE CURRENT passes, in cm, and

U_{peak} = peak HF test voltage.

The corresponding limit for cords intended for BIPOLAR application is

$$I_{\text{leakage}} [\text{mA}] = 1,8 \times 10^{-5} \times d \times L \times f_{\text{test}} \times U_{\text{peak}}$$

Compliance is checked as follows:

The full length of the sample insulation except that within 1 cm of exposed conductors, but no more than 30 cm of length, shall be immersed in a 0,9 % saline bath or wrapped in a saline-soaked porous cloth. All operative inner conductors shall be connected together to one pole of an HF voltage source having an approximately sinusoidal waveform and a frequency f_{test} of 300 kHz to 1 MHz. The opposite pole of the HF voltage source is connected to a conductive electrode immersed in the saline bath or to foil wrapped around the midsection of the saline-soaked cloth. HF LEAKAGE CURRENT I_{leakage} is monitored by means of a suitable instrument connected in series with the HF voltage source output. The HF test voltage U_{peak} is monitored between the HF voltage source output poles.

The HF test voltage U_{peak} is advanced until the peak voltage equals the lesser of RATED ACCESSORY VOLTAGE or 400 V_{peak}. The measured HF LEAKAGE CURRENT I_{leakage} shall not exceed the specified limit.

b) Measured HF leakage capacitance

The preceding item a) may alternatively be fulfilled by limiting measured HF leakage capacitance for cords intended for MONOPOLAR application to no more than

$$C_{\text{leakage}} [\text{pF}] = 2 \times d \times L$$

and for cords intended for BIPOLAR application to no more than

$$C_{\text{leakage}} [\text{pF}] = 4 \times d \times L$$

where:

d = smallest outer dimension of the insulation, in mm, and

L = length of sample insulation immersed in saline bath, in cm.

The measured HF leakage capacitance shall not exceed the specified relevant limit.

Customer: Bradley Hyde- No. of User(s): 1 - Company: Latham and Watkins

Order No.: WS-2023-008299 - IMPORTANT: This file is copyright of IEC, Geneva, Switzerland. All rights reserved.

This file is subject to a licence agreement. Enquiries to Email: sales@iec.ch - Tel.: +41 22 919 02 11

Compliance is checked as follows:

The full length of the sample insulation except that within 1 cm of exposed conductors, but no more than 30 cm of length, shall be immersed in a 0,9 % saline bath. All operative inner conductors shall be connected together to one measuring terminal of a capacitance-measuring instrument having a sensing frequency of 100 kHz to 1 MHz. The opposite measuring terminal of the capacitance measuring instrument is connected to a conductive electrode immersed in the saline bath. HF leakage capacitance is the capacitance indicated by the capacitance measuring instrument when operated according the instrument manufacturer's recommended practices.

NOTE HF leakage limits and tests for all parts of ACTIVE ACCESSORIES are under consideration.

201.8.8.3.103 * ACTIVE ACCESSORY HF dielectric strength

The insulation applied to ACTIVE ACCESSORIES shall be capable of withstanding HF voltage of 120 % of the RATED ACCESSORY VOLTAGE.

Compliance is checked as follows:

The tests shall be performed at a test voltage related to the RATED ACCESSORY VOLTAGE specified by the manufacturer of the HF SURGICAL ACCESSORY in the instructions for use (see 201.7.9.2.14 e), as detailed in the following test methods. For the cords of ACTIVE ACCESSORIES and ACTIVE ELECTRODES, a portion of the insulation which has been preconditioned in saline is wound with a maximum of five turns of bare conductive wire having a diameter of 0,4 mm $\pm$ 10 % at a pitch of at least 3 mm without deforming the surface of the sample. If necessary to prevent inadvertent arc discharge, the CREEPAGE DISTANCE between this wire and operative conductive parts of ACTIVE ELECTRODES may be increased to 10 mm by application of insulation. Such added insulation shall have a thickness no greater than 1 mm and shall cover no more than 2 mm of ACTIVE ELECTRODE INSULATION. One pole of the HF test voltage source shall be connected to the bare conductive test wire, and the opposite pole shall be connected simultaneously to all operative conductors in the sample being tested.

ACTIVE HANDLES, together with any detachable cords and detachable ACTIVE ELECTRODES which are specified as compatible, shall be wrapped in a porous cloth soaked in 0,9 % saline. This cloth shall cover the entire exterior surface of the handle and extend at least 150 mm on to the surface of the cord and 5 mm on to the ACTIVE ELECTRODE INSULATION. If necessary, the CREEPAGE DISTANCE between the cloth and exposed operative conductive parts of the ACTIVE ELECTRODE may be insulated as described above. The midsection of the saline-soaked cloth is wrapped with metal foil and connected to one pole of the HF test voltage source. All operative inner conductors in the samples being tested, including the operative tip(s) of the ACTIVE ELECTRODE, shall be connected simultaneously to the opposite pole.

The peak HF test voltage is monitored between the HF voltage source output poles. The output of the HF test voltage source is then increased until the peak voltage equals 120 % of the peak voltage according to the RATED ACCESSORY VOLTAGE and maintained for 30 s in such a manner that it stresses the insulation of the test sample. No breakdown of the insulation material shall occur and the same insulation shall subsequently be tested at mains frequency according to 201.8.8.3.104.

NOTE Blue corona is normal and is not considered a breakdown of insulation.

Those parts of the test samples which are not insulated in NORMAL USE shall be adequately protected against contact with the saline solution during preconditioning, and this protection shall be left in place during the tests.

Tests:

Apply an approximately sinusoidal voltage at a frequency of 400 kHz $\pm$ 100 kHz with a continuous waveform, or alternately with a modulated waveform (modulation frequency higher

Customer: Bradley Hyde- No. of User(s): 1 - Company: Latham and Watkins

Order No.: WS-2023-008299 - IMPORTANT: This file is copyright of IEC, Geneva, Switzerland. All rights reserved.

This file is subject to a licence agreement. Enquiries to Email: sales@iec.ch - Tel.: +41 22 919 02 11

than 10 kHz) with the peak test voltage equal to 120 % of the peak voltage according to the RATED ACCESSORY VOLTAGE specified by the manufacturer of the HF SURGICAL ACCESSORY and with a test CREST FACTOR (cf_{test}) which is defined as follows:

For RATED ACCESSORY VOLTAGES less than or equal 1 600 V:

$$cf_{test} \leq 2$$

For RATED ACCESSORY VOLTAGES greater 1 600 V and less than or equal to 4 000 V:

$$cf_{test} = \frac{U_{acc} - 400 [\text{Volts}]}{600 [\text{Volts}]} \pm 10 \% \text{ where}$$

U_{acc} = RATED ACCESSORY VOLTAGE in volts

For RATED ACCESSORY VOLTAGES greater 4 000 V:

$$cf_{test} = 6 \pm 10 \%$$

ACTIVE ACCESSORIES intended to be used with HF SURGICAL MODES or output settings requiring specific approval shall withstand 120 % of the peak output voltage of such HF SURGICAL MODE or output setting. They shall be tested under the same conditions as described above but with the actual crest factor of such HF SURGICAL MODE or output setting (see 201.7.9.2.2.101 c) 3)).

201.8.8.3.104 * ACTIVE ACCESSORY mains frequency dielectric strength

The insulation applied to an ACTIVE ACCESSORY, including those portions of insulation having been tested at HF according to 201.8.8.3.103, shall withstand a d.c. or mains frequency peak voltage of 1 000 V greater than the RATED ACCESSORY VOLTAGE specified by the manufacturer of the HF SURGICAL ACCESSORY.

Compliance is tested as follows:

The test voltage source shall produce a d.c. or mains frequency signal. The test duration shall be 30 s for ACTIVE HANDLES and ACTIVE CONNECTORS. The test duration for the cords of ACTIVE ACCESSORIES shall be 5 min. Although corona discharge may occur, no breakdown of the insulation or flashover shall occur. Immediately after this dielectric strength test, any incorporated FINGERSWITCH shall be operated 10 times. An ohmmeter shall be used to test whether the switching mechanism operates as intended to ensure that, when connected to HF SURGICAL EQUIPMENT, the HF output will be de-energized when the FINGERSWITCH is released.

The insulated parts of ACTIVE CONNECTORS more than 10 mm CREEPAGE DISTANCE from exposed operative conductors shall be wrapped with a porous cloth soaked in 0,9 % saline. The midsection of the cloth is then wrapped with metal foil. The test voltage is applied between the foil and all of the operative ACTIVE CONNECTOR contacts.

The entire length of the insulation of cords of ACTIVE ACCESSORIES, including that portion previously tested at HF according to 201.8.8.3.103, but exclusive of the sections within 100 mm of the ends, shall be immersed in a bath of 0,9 % saline. The test voltage is applied between a conductive electrode immersed in the saline bath and all of the conductors in the cord simultaneously.

ACTIVE HANDLES complete with detachable electrodes are prepared for testing and connected to the test voltage source using the same techniques as described in 201.8.8.3.103. The saline-soaked cloth and foil applied for that test may be left in place for this test provided care is taken to ensure that the cloth remains thoroughly wetted.

Customer: Bradley Hyde- No. of User(s): 1 - Company: Latham and Watkins

Order No.: WS-2023-008299 - IMPORTANT: This file is copyright of IEC, Geneva, Switzerland. All rights reserved.

This file is subject to a licence agreement. Enquiries to Email: sales@iec.ch - Tel.: +41 22 919 02 11

201.8.9.1.5 ME EQUIPMENT RATED for high altitudes

Amendment:

This requirement does not apply for the separation between HF PATIENT CIRCUITS and the ENCLOSURE including SIGNAL INPUT PARTS and SIGNAL OUTPUT PARTS and between different HF PATIENT CIRCUITS.

For HF SURGICAL EQUIPMENT and ASSOCIATED EQUIPMENT, requirements for separation between HF PATIENT CIRCUITS and the ENCLOSURE including SIGNAL INPUT PARTS and SIGNAL OUTPUT PARTS and between different HF PATIENT CIRCUITS are specified by 201.8.5.1.2.

201.8.10.4 Cord-connected HAND-HELD parts and cord-connected foot-operated control devices**201.8.10.4.1 Limitation of operating voltages**

Subclause 8.10.4.1 of the general standard does not apply. See 201.8.10.4.101.

201.8.10.4.2 * Connection cords

Replacement:

Anchorage of cords of ACTIVE ACCESSORIES shall be designed to minimize the risk to PATIENTS and OPERATORS arising from damage to conductors or insulation caused by cable flexure or excessive tension.

Compliance shall be checked by inspection and by the following test:

The anchorages on ACTIVE HANDLES and ACTIVE CONNECTORS are tested one at a time.

The ACTIVE HANDLE or ACTIVE CONNECTOR under test is FIXED in an apparatus similar to that shown in Figure 201.108, so that when the oscillating member of the apparatus is at the middle of its travel, the axis of the cord, where it leaves the part under test, is vertical and passes through the axis of oscillation. The cord is passed through an aperture 300 mm from the axis of oscillation and a weight equal to the cord and connector of the ACTIVE ACCESSORY is affixed to the cable below this aperture for the purpose of applying tension to the cord. The maximum diameter of the hole should not be more than 2 times the diameter of the cord.

Where an anchorage of the ACTIVE HANDLE or ACTIVE CONNECTOR under test is fitted with two or more cords, these shall be tested together, with the total weight affixed to the anchorage being the sum of the weights required to be applied to each cord individually.

The oscillating member is rotated through an angle of 90 ° (45 ° on each side of the vertical).

The number of cycles applied to cable anchorages of ACTIVE HANDLES shall be 10 000 (200 for ACTIVE ACCESSORIES marked for single use only) at the rate of approximately 30 cycles per minute. The number of cycles applied to anchorages of cables of ACTIVE CONNECTORS shall be 5 000 (100 for ACTIVE ACCESSORIES marked for single use only) at the rate of approximately 30 cycles per minute.

After the test, the cord shall not have worked loose nor shall it show any damage. For multi-conductor cables there shall be no short circuits between individual conductors. The tensioning weight shall be increased to 1 kg and individual conductors checked for continuity using a d.c., current not in excess of 1 A.

Customer: Bradley Hyde- No. of User(s): 1 - Company: Latham and Watkins

Order No.: WS-2023-008299 - IMPORTANT: This file is copyright of IEC, Geneva, Switzerland. All rights reserved.

This file is subject to a licence agreement. Enquiries to Email: sales@iec.ch - Tel.: +41 22 919 02 11

SERENDIA-00084181

Ex. 2077-033

CX-2848.33

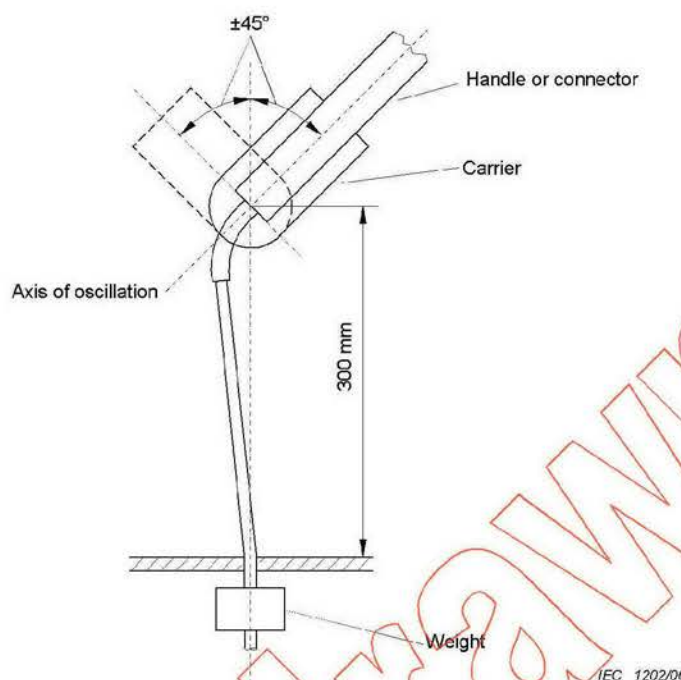


Figure 201.108 – Test apparatus for anchorages of cords of ACTIVE ACCESSORY

Additional subclauses:

201.8.10.4.101 * SWITCH SENSORS

201.8.10.4.101.1 General

Except where provided for in subclause 201.8.10.4.101.2, HF SURGICAL EQUIPMENT and applicable ASSOCIATED EQUIPMENT shall be provided with a SWITCH SENSOR requiring continuous activation in order to energize the ACTIVE OUTPUT TERMINALS.

The SWITCH SENSOR shall be supplied from a power source isolated from the MAINS PART and from earth, having a voltage not exceeding 12 V, if a CONDUCTIVE CONNECTION to the APPLIED PART exists, and not exceeding 24 V a.c. or 34 V d.c. in other cases.

NOTE This requirement applies to voltages appearing within SWITCH SENSORS. Common-mode HF voltages should be disregarded.

Under SINGLE FAULT CONDITION the SWITCH SENSOR shall not cause low-frequency PATIENT LEAKAGE CURRENT(S) exceeding the allowable limits (see 201.8.7.3).

Compliance is checked by inspection, functional check, and by measurement of voltage and LEAKAGE CURRENT(S).

Where the SWITCH SENSOR is provided with input terminals intended for connection to external electrical switch contacts, it shall not be possible to activate any output of the HF SURGICAL EQUIPMENT when the input terminals are bridged by a resistance equal to or greater than 1 000 Ω .

Compliance is checked by a functional test.

Customer: Bradley Hyde- No. of User(s): 1 - Company: Latham and Watkins

Order No.: WS-2023-008299 - IMPORTANT: This file is copyright of IEC, Geneva, Switzerland. All rights reserved.

This file is subject to a licence agreement. Enquiries to Email: sales@iec.ch - Tel.: +41 22 919 02 11

Each SWITCH SENSOR shall activate only its intended single ACTIVE OUTPUT TERMINAL and shall control no more than one HF SURGICAL MODE at any one time.

NOTE For the purpose of this requirement the two arms of a rocker style switch are considered to be two individual switches.

201.8.10.4.101.2 Non-continuous activation

Non-continuous activation mode of the SWITCH SENSOR is accepted only if

- a) the output of the HF SURGICAL EQUIPMENT is automatically stopped in accordance with the specific application of the equipment;
- b) a visible indicator is provided to indicate to the OPERATOR that the HF SURGICAL EQUIPMENT is set to such a specific application mode, and
- c) a means of manual output deactivation is provided.

Compliance is checked by inspection of ACCOMPANYING DOCUMENTS and functional test.

201.8.10.4.101.3 Impedance sensing activation

A SWITCH SENSOR which is intended to activate HF output in response to the impedance appearing between BIPOLAR ACTIVE OUTPUT TERMINALS is accepted only for BIPOLAR COAGULATION.

Where such an impedance-sensing SWITCH SENSOR is provided as an alternative or in addition to a contact-closure sensing SWITCH SENSOR, then

- a) it shall not be possible under any conditions for HF output to be energized solely as a result of interruption and restoration of the SUPPLY MAINS, and
- b) impedance-sensing activation shall be enabled only in response to a specific OPERATOR selection, and
- c) that selection shall be visibly indicated to the OPERATOR.

Impedance sensing SWITCH SENSORS shall not be permitted for MONOPOLAR HF output activation. The requirements of this subclause do not apply to SWITCH SENSORS which are capable only of automatically terminating HF output according to the purpose of specific application modes (see 201.8.10.4.101.2 a)).

Compliance is checked by inspection of ACCOMPANYING DOCUMENTS and functional test.

201.8.10.4.101.4 Footswitches

Footswitches shall comply with the following requirement (see also 201.11.6.5 and 201.12.2).

The force required to actuate the switch shall be not less than 10 N, applied over an area of 625 mm² anywhere on the operating surface of the footswitch.

Compliance is checked by measurement of the actuating force.

201.9 Protection against mechanical hazards of ME EQUIPMENT and ME SYSTEMS

Clause 9 of the general standard applies.

201.10 Protection against unwanted and excessive radiation HAZARDS

Clause 10 of the general standard applies.

Customer: Bradley Hyde- No. of User(s): 1 - Company: Latham and Watkins

Order No.: WS-2023-008299 - IMPORTANT: This file is copyright of IEC, Geneva, Switzerland. All rights reserved.

This file is subject to a licence agreement. Enquiries to Email: sales@iec.ch - Tel.: +41 22 919 02 11

201.11 Protection against excessive temperatures and other HAZARDS

Clause 11 of the general standard applies, except as follows:

201.11.1.1 * Maximum temperature during NORMAL USE

Addition

Duty cycle: HF SURGICAL EQUIPMENT, set up to deliver its RATED OUTPUT POWER into a resistive load using the electrode cable, is operated for 1 h with a DUTY CYCLE as specified by the manufacturer but with operating times of at least 10 s alternating with a resting time of not more than 30 s.

201.11.1.2.1 APPLIED PARTS intended to supply heat to a PATIENT

Addition:

ACTIVE ELECTRODES are considered to be APPLIED PARTS intended to supply heat to a PATIENT as part of their intended clinical effect (CUTTING and COAGULATION). Disclosure of temperatures and clinical effects is not required.

201.11.1.2.2 APPLIED PARTS not intended to supply heat to a PATIENT

Addition:

NEUTRAL ELECTRODES are considered to be APPLIED PARTS not intended to supply heat to a PATIENT.

201.11.6.3 * Spillage on ME EQUIPMENT and ME SYSTEM

Replacement:

The ENCLOSURE of the HF SURGICAL EQUIPMENT and ASSOCIATED EQUIPMENT shall be constructed so that liquid spillage in NORMAL USE does not wet electrical insulation or other components which, when wetted, are likely to affect adversely the safety of the HF SURGICAL EQUIPMENT and ASSOCIATED EQUIPMENT.

Compliance is checked by the following test.

A quantity of one litre of water is poured steadily onto the middle of the top surface of the HF SURGICAL EQUIPMENT and ASSOCIATED EQUIPMENT over a period of 15 s. HF SURGICAL EQUIPMENT and ASSOCIATED EQUIPMENT intended to be built into a wall or cabinet is tested mounted as recommended, the water being poured onto the wall above the control panel. After this treatment, the HF SURGICAL EQUIPMENT and ASSOCIATED EQUIPMENT shall withstand the dielectric strength test specified in 201.8.8.3, and inspection shall show that water which may have entered the ENCLOSURE cannot adversely affect the safety of the HF SURGICAL EQUIPMENT and ASSOCIATED EQUIPMENT. In particular, there shall be no trace of water on the insulation for which CREEPAGE DISTANCES are specified in 8.9.1 of the general standard.

201.11.6.5 Ingress of water or particulate matter into ME EQUIPMENT and ME SYSTEMS

Addition:

- a) * The electrical switching parts of footswitches for HF SURGICAL EQUIPMENT and ASSOCIATED EQUIPMENT intended for use in operating rooms shall be protected against the effects of ingress of liquids that might cause inadvertent energization of the APPLIED PART.

Compliance is checked by the following test.

Customer: Bradley Hyde- No. of User(s): 1 - Company: Latham and Watkins

Order No.: WS-2023-008299 - IMPORTANT: This file is copyright of IEC, Geneva, Switzerland. All rights reserved.

This file is subject to a licence agreement. Enquiries to Email: sales@iec.ch - Tel.: +41 22 919 02 11

The footswitch shall be completely immersed in 0,9 % saline to a depth of 150 mm for a period of 30 min. While immersed, it shall be connected to a SWITCH SENSOR corresponding to its NORMAL USE and actuated 50 times. The SWITCH SENSOR shall register deactivation upon each release.

- b) * The electrical parts of FINGERSWITCHES shall be protected against the effects of ingress of liquids that might cause inadvertent energization of the APPLIED PART (see also 201.8.8.3.103).

Compliance is checked by the following test.

The a.c. impedance of each of the switching terminals of the ACTIVE CONNECTOR shall be measured using a frequency of at least 1 kHz and a voltage of less than 12 V. The ACTIVE HANDLE is supported horizontally at least 50 mm above any surface with the switch activating parts uppermost. One litre of 0,9 % saline solution is poured steadily from above over the ACTIVE HANDLE over a period of 15 s so as to wet the entire length of the ACTIVE HANDLE. The liquid is allowed to drain away freely. The a.c. impedance of the switching terminals shall remain greater than 2 000 Ω .

Immediately after, each FINGERSWITCH is operated and released 10 times. The a.c. impedance of the switching terminals shall exceed 2 000 Ω within 0,5 s after each release.

201.11.6.7 * Sterilization of ME EQUIPMENT and ME SYSTEMS

Addition:

Unless marked for single use only, ACTIVE ACCESSORIES and all detachable parts thereof, except ACTIVE CONNECTORS detachable from cords without use of tools, shall comply with the requirements of this particular standard after being tested according to this subclause of the general standard.

201.11.8 Interruption of the power supply / SUPPLY MAINS to ME EQUIPMENT

Addition:

When HF SURGICAL EQUIPMENT is switched off and on again or when the SUPPLY MAINS is interrupted and re-established

- the output power for a given setting of the output control shall not increase by more than 20 %, and
- the HF SURGICAL MODE shall not be changed except to a stand-by mode in which no output is produced.

Compliance is checked by measurement of the power, averaged over a period of 1 s, and observation of the operating mode

- a) *with repeated operation of the mains switch of the HF SURGICAL EQUIPMENT;*
- b) *with interruption and re-establishment of the SUPPLY MAINS, the switch in the HF SURGICAL EQUIPMENT being left in the "ON" position.*

201.12 Accuracy of controls and instruments and protection against hazardous outputs

Clause 12 of the general standard applies, except as follows:

201.12.1 Accuracy of controls and instruments

Additional subclauses:

Customer: Bradley Hyde- No. of User(s): 1 - Company: Latham and Watkins

Order No.: WS-2023-008299 - IMPORTANT: This file is copyright of IEC, Geneva, Switzerland. All rights reserved.

This file is subject to a licence agreement. Enquiries to Email: sales@iec.ch - Tel.: +41 22 919 02 11

201.12.1.101 Accuracy of output control setting

For output powers in excess of 10 % of the RATED OUTPUT POWER, the actual power as a function of the load resistance and output control setting shall not deviate from that shown in the diagrams specified in 201.7.9.3.1 by more than ± 20 %.

Compliance is checked by performing the test of 201.12.1.102 but using appropriate values of load resistance.


201.12.1.102 Monotonicity of output control setting

The output power shall not increase with the decrease of the output control setting (see 201.7.9.3.1, Figure 201.109 and Figure 201.110).

Compliance is checked by the following test:

a) * MONOPOLAR outputs

The output power as a function of the output control setting is measured at a minimum of five particular values of the load resistance, including 100 Ω , 200 Ω , 500 Ω , 1 000 Ω , 2 000 Ω and at the RATED LOAD. ACTIVE ACCESSORIES and NEUTRAL ELECTRODES supplied with HF SURGICAL EQUIPMENT or 3 m lengths of insulated conductors shall be used for connection of the load resistors.

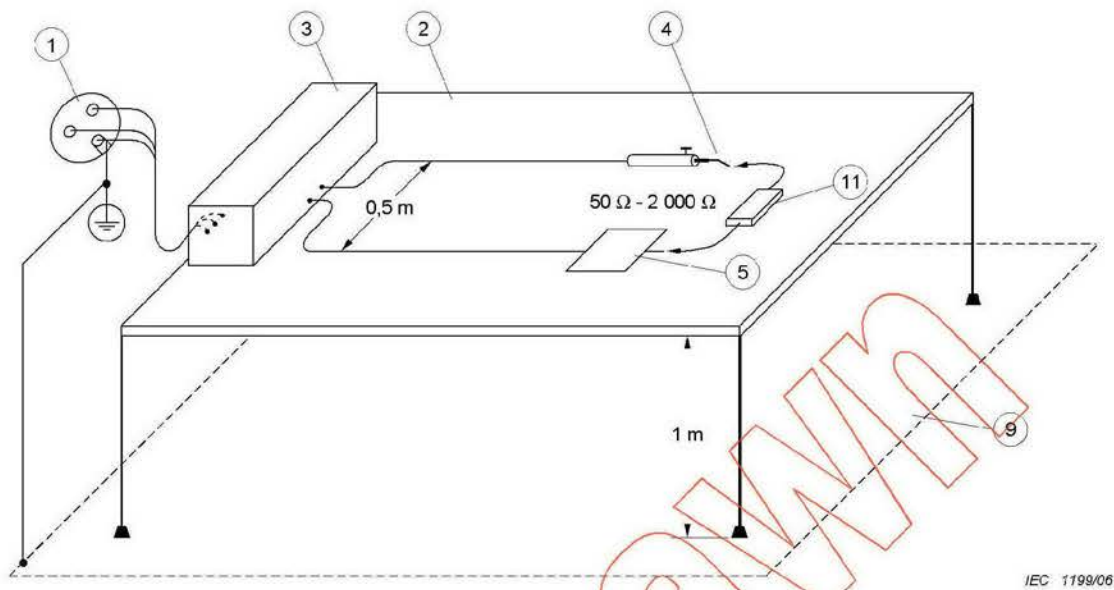

IEC 1199/06

Key


- ① SUPPLY MAINS
- ② Table, made of insulating material
- ③ HF SURGICAL EQUIPMENT
- ④ ACTIVE ELECTRODE
- ⑤ NEUTRAL ELECTRODE, metallic or in contact with metal foil of the same size
- ⑨ Earthed conductive plane
- ⑪ Load resistance as required with HF power measuring device

Figure 201.109 – Measurement of output power – MONOPOLAR output

b) * BIPOLAR outputs

The output power as a function of the output control setting is measured at a minimum of five particular values of the load resistance, including 10 Ω , 50 Ω , 200 Ω , 500 Ω , 1 000 Ω and at the RATED LOAD. The BIPOLAR cord supplied with the HF SURGICAL EQUIPMENT or a 3 m length of two conductor insulated cord RATED 600 V or greater shall be used for the connection of the load resistors.

MANUFACTURERS shall provide specific instructions on how to set up these measurements on alternate forms of BIPOLAR ELECTRODES.

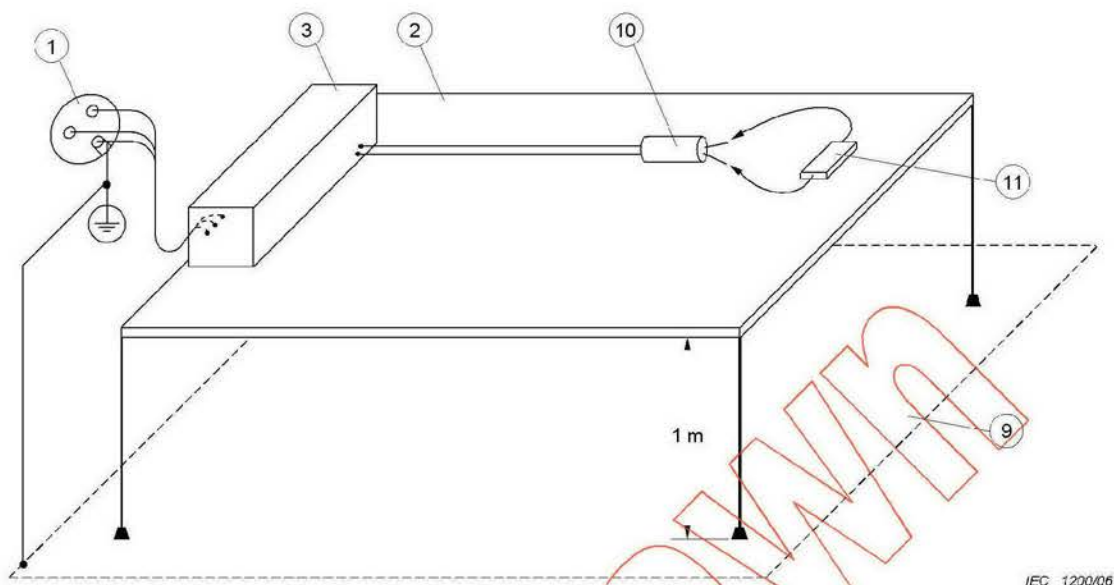

IEC 1200/06

Key

- ① SUPPLY MAINS
- ② Table, made of insulating material
- ③ HF SURGICAL EQUIPMENT
- ⑨ Earthed conductive plane
- ⑩ Activated BIPOLAR ELECTRODE
- ⑪ Load resistance as required with HF power measuring device

Figure 201.110 – Measurement of output power – BIPOLAR output

201.12.1.103 * Accuracy of MAXIMUM OUTPUT VOLTAGE

For each HF SURGICAL MODE available in HF SURGICAL EQUIPMENT, the MAXIMUM OUTPUT VOLTAGE applied to the ACTIVE OUTPUT TERMINALS shall not exceed that specified in 201.7.9.3.1.

Compliance is checked by use of an oscilloscope. See also 201.5.4 aa). Measurements shall be taken at the output setting and load condition which yields the highest peak output voltage for each HF SURGICAL MODE.

201.12.2 Usability

Addition:

- a) Where a double footswitch assembly is used to select CUTTING and COAGULATION output modes, the arrangement shall be such that, when viewed by the OPERATOR, the left pedal activates CUTTING and the right pedal activates COAGULATION.

Compliance is checked by inspection.

- b) * In an ACTIVE HANDLE which incorporates separate FINGERSWITCHES for selectively activating CUTTING and COAGULATION HF SURGICAL MODES, that which activates CUTTING shall be nearer to the ACTIVE ELECTRODE than is the other.

Compliance is checked by inspection.

- c) It shall not be possible to energize simultaneously more than one ACTIVE OUTPUT TERMINAL, unless:

Customer: Bradley Hyde- No. of User(s): 1 - Company: Latham and Watkins

Order No.: WS-2023-008299 - IMPORTANT: This file is copyright of IEC, Geneva, Switzerland. All rights reserved.

This file is subject to a licence agreement. Enquiries to Email: sales@iec.ch - Tel.: +41 22 919 02 11

SERENDIA-00084188

Ex. 2077-040

CX-2848.40

- 1) each ACTIVE OUTPUT TERMINAL has independent sets of controls for selection of HF SURGICAL MODE, HF output setting and independent SWITCH SENSORS,

or

- 2) two MONOPOLAR ACTIVE OUTPUT TERMINALS have independent SWITCH SENSORS and share a common FULGURATION output.

Compliance is checked by inspection and functional check.

- d) * During simultaneous activation the audible tone shall be different from the tone produced during single output activation. See also 201.12.4.2.101. Under no circumstances shall any PATIENT circuit become energized by more than is defined in 201.8.7.3.101 c), unless the output for that PATIENT circuit is activated by the OPERATOR.

Compliance is checked by inspection and functional check.

- e) * ACTIVE OUTPUT TERMINALS on HF SURGICAL EQUIPMENT and ASSOCIATED EQUIPMENT shall differ in configuration sufficiently such that MONOPOLAR ACTIVE ACCESSORIES, NEUTRAL ELECTRODES and BIPOLAR ACTIVE ACCESSORIES cannot be improperly connected.

NOTE See Annex AA.

Compliance is checked by inspection.

- f) * ACTIVE CONNECTORS having more than one pin shall have FIXED pin spacing. "Flying leads" are prohibited.

Compliance is checked by inspection.

- g) * Where more than one HF SURGICAL MODE can be energized by a single SWITCH SENSOR, an indication shall be provided to show which HF SURGICAL MODE is selected before an output is energized.

Compliance is checked by inspection and functional test.

201.12.4.2 * Indication of parameters relevant to safety

Addition:

If the total output power in any HF SURGICAL MODE, including simultaneous activation of independent outputs if available, exceeds 400 W averaged over any period of 1 s when each of the outputs is terminated at the RATED LOAD, then special consideration of potential HAZARDS shall be addressed in the RISK MANAGEMENT FILE, especially with regard to NEUTRAL ELECTRODES.

Compliance is checked by measurement.

Additional subclause:

201.12.4.2.101 Output indicator

HF SURGICAL EQUIPMENT shall be provided with a device which gives an audible signal when any output circuit is energized by the operation of a SWITCH SENSOR or as a result of a SINGLE FAULT CONDITION. The sound output shall have its major energy content in the band of frequencies between 100 Hz and 3 kHz. The sound source shall be capable of producing a sound level of at least 65 dBA at a distance of 1 m from the HF SURGICAL EQUIPMENT according to the one direction specified by the manufacturer. An accessible sound level control may be provided, but shall not reduce the sound level below 40 dBA. For simultaneous activation see also 201.12.2 d).

In order that the OPERATOR may distinguish between the audible alarm called for in 201.8.4.101 and the signal specified above, either the former shall be pulsed or two different frequencies shall be employed.

NOTE This audible signal is not intended to meet the definition of ALARM SIGNAL in IEC 60601-1-8. See also Clause 208 of this standard.

Customer: Bradley Hyde- No. of User(s): 1 - Company: Latham and Watkins

Order No.: WS-2023-008299 - IMPORTANT: This file is copyright of IEC, Geneva, Switzerland. All rights reserved.

This file is subject to a licence agreement. Enquiries to Email: sales@iec.ch - Tel.: +41 22 919 02 11

Compliance is checked by functional check and measurement of the sound level.

201.12.4.3 Accidental selection of excessive output values

Additional subclause:

201.12.4.3.101 * Output reduction means

Except as provided for in 201.7.9.2.2.101 a) item 7 and 201.7.9.3.1, for each HF SURGICAL MODE, HF SURGICAL EQUIPMENT shall incorporate means to enable the output power to be reduced to not more than 5 % of the RATED OUTPUT POWER or 10 W, whichever is smaller (see also 201.12.1.102).

Compliance is checked by measurement of output power and inspection.

201.12.4.4 Incorrect output

Additional subclauses:

201.12.4.4.101 * Maximum allowed output power in SINGLE FAULT CONDITIONS

MONOPOLAR HF SURGICAL EQUIPMENT having a RATED OUTPUT POWER greater than 50 W and all BIPOLAR outputs of HF EQUIPMENT shall be provided with an alarm and/or interlock system to indicate and/or prevent a significant increase in the output power relative to the output setting.

The maximum allowed output power under SINGLE FAULT CONDITIONS shall be calculated separately for each PATIENT CIRCUIT and operation mode.

The maximum allowed output power in SINGLE FAULT CONDITIONS is defined according to Table 201.102:

Table 201.102 – Maximum output powers in SINGLE FAULT CONDITIONS

Setting (range in % of RATED OUTPUT POWER)	Maximum allowed output power in SINGLE FAULT CONDITIONS
Less than 10	20 % of RATED OUTPUT POWER
10 to 25	Setting x 2
Greater than 25 and up to 80	Setting + 25 % of RATED OUTPUT POWER
Greater than 80 and up to 100	Setting + 30 % of RATED OUTPUT POWER

Compliance is checked by examination of the technical documentation and testing by simulation of appropriate SINGLE FAULT CONDITIONS.

201.12.4.4.102 * Output power during simultaneous activation

For HF SURGICAL EQUIPMENT providing simultaneous activation of more than one PATIENT circuit (see 201.12.2), the PATIENT circuits shall not deliver an output power that exceeds the range of deviation defined in 201.12.1.101 by more than 20 % when they are simultaneously activated under any available combination of HF SURGICAL MODES.

Any single activated PATIENT circuit shall comply with 201.12.1.101.

Compliance is checked by the following tests (see Figure 201.111).

For HF SURGICAL EQUIPMENT as defined in 201.12.2 c):

Customer: Bradley Hyde- No. of User(s): 1 - Company: Latham and Watkins
Order No.: WS-2023-008299 - IMPORTANT: This file is copyright of IEC, Geneva, Switzerland. All rights reserved.
This file is subject to a licence agreement. Enquiries to Email: sales@iec.ch - Tel.: +41 22 919 02 11

The output under test is activated at 20 % of its RATED OUTPUT POWER and the HF current reading of this output noted. Any other output is then activated at maximum power and the current of the output under test shall not increase by more than 10 %.

The output under test is activated at 50 % and at 100 % output settings and the current values noted. These values shall not increase by more than 10 % when the other output is activated additionally.

These tests are repeated with all possible combinations of outputs which may be activated together at any one time.

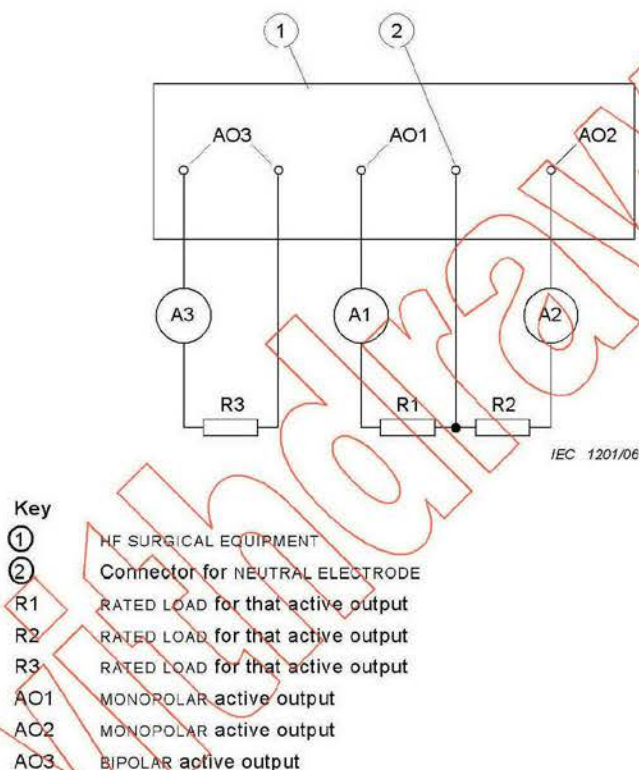


Figure 201.111 – Method of testing feedback from one active output to another in simultaneous activation

201.13 HAZARDOUS SITUATIONS and fault conditions

Clause 13 of the general standard applies, except as follows:

201.13.2.13 Overload

Additional subclause:

201.13.2.13.101 * Protection against the effects of short-circuiting of the electrodes

HF SURGICAL EQUIPMENT shall be capable of withstanding, without damage, the effects of short-circuiting or open-circuiting the output when energized at maximum output setting.

Compliance is checked by the following test.

Customer: Bradley Hyde- No. of User(s): 1 - Company: Latham and Watkins

Order No.: WS-2023-008299 - IMPORTANT: This file is copyright of IEC, Geneva, Switzerland. All rights reserved.

This file is subject to a licence agreement. Enquiries to Email: sales@iec.ch - Tel.: +41 22 919 02 11

Connect the conductors described in 201.12.1.102, items a) and b), to the PATIENT circuit connections and, for each HF SURGICAL MODE, set the output control to the maximum position. The output is then switched on, and the remote ends of the activated pair of conductors are short-circuited for a period of 5 s and then open-circuited for a period of 15 s. The output is then switched off for a period of 1 min. The above cycle is repeated for a total of 10 times.

After this test the HF SURGICAL EQUIPMENT shall comply with all the requirements of this particular standard.

201.14 PROGRAMMABLE ELECTRICAL MEDICAL SYSTEMS (PEMS)

Clause 14 of the general standard applies.

201.15 Construction of ME EQUIPMENT

Clause 15 of the general standard applies, except as follows:

201.15.4.1 Construction of connectors

Additional subclauses:

201.15.4.1.101 * Compatibility with third party ACTIVE ELECTRODES

The manufacturer of an ACTIVE ACCESSORY with a detachable ACTIVE ELECTRODE shall provide upon request the dimensions and associated tolerances for the mating part of any ACTIVE ELECTRODE which is intended to be attached to the ACTIVE ACCESSORY.

Compliance is checked by inspection of the ACCOMPANYING DOCUMENTS.

The manufacturer of an ACTIVE ACCESSORY with a detachable ACTIVE ELECTRODE shall specify in the ACCOMPANYING DOCUMENTS the ACTIVE ELECTRODES with which it is intended to be compatible.

Compliance is checked by demonstrating conformance with all relevant requirements of this particular standard.

201.15.4.1.102 * Retention of detachable ACTIVE ELECTRODES

The manufacturer of a detachable ACTIVE ELECTRODE shall specify in its ACCOMPANYING DOCUMENTS the ACTIVE ACCESSORIES with which it is intended to be used.

The detachable ACTIVE ELECTRODE shall fit securely into the specified ACTIVE ACCESSORIES.

Compliance shall be checked by inspection and by the following test:

The detachable ACTIVE ELECTRODE is inserted ten times into a specified ACTIVE ACCESSORY. Afterwards, the ACTIVE ELECTRODE shall not detach when subjected to a pull equivalent to ten times the weight of the ACTIVE ELECTRODE up to a maximum of 10 N for one minute along the axis of insertion.

When a detachable ACTIVE ELECTRODE is inserted into a specified ACTIVE ACCESSORY, the combination shall conform to all other applicable requirements of this particular standard.

Additional subclauses:

Customer: Bradley Hyde- No. of User(s): 1 - Company: Latham and Watkins

Order No.: WS-2023-008299 - IMPORTANT: This file is copyright of IEC, Geneva, Switzerland. All rights reserved.

This file is subject to a licence agreement. Enquiries to Email: sales@iec.ch - Tel.: +41 22 919 02 11

201.15.101 * NEUTRAL ELECTRODES**201.15.101.1 General requirements for NEUTRAL ELECTRODES**

Except for any PATIENT circuit intended only for connection to a BIPOLAR ELECTRODE, HF SURGICAL EQUIPMENT having a RATED OUTPUT POWER in excess of 50 W shall be provided with a NEUTRAL ELECTRODE.

Compliance is checked by inspection.

201.15.101.2 * NE cord attachment

The NEUTRAL ELECTRODE shall be reliably connected to the cord. Except for a MONITORING NE, any current used for monitoring the electrical continuity of the electrode cord and its connections shall pass through a section of the electrode.

Compliance is checked by the following test.

An electrical continuity test is conducted using a current of at least 1 A but not more than 5 A from a d.c. or mains frequency current source with a no-load voltage not exceeding 6 V. The resistance shall be 1 Ω or less.

201.15.101.3 * NE cord connector, no conductive parts on PATIENT

Any contacts of the electrical connector of an NE cord for attachment to a detachable NE shall be designed so that their conductive parts cannot come into contact with the body of the PATIENT in the event of inadvertent disconnection.

Compliance is checked by the following test.

The NE cord is detached from the NE and, using the standard test finger shown in Figure 6 of the general standard, it is verified that contact with conductive parts of the cable connector is not possible.

201.15.101.4 * NE cord insulation

The insulation of NE cords shall be adequate to prevent a burn injury to the PATIENT and the OPERATOR.

Compliance is checked by application of the following tests in the order shown:

- *HF leakage test according to 201.8.8.3.102 a) with a test voltage [U_{peak}] of 400 V_{peak} . HF LEAKAGE CURRENT shall not exceed*

$$I_{leakage} [mA] = 1,8 \times 10^{-5} \times d \times L \times f_{test} \times U_{peak}$$

or alternatively the HF leakage capacitance test according to 201.8.8.3.102 b). The HF leakage capacitance shall not exceed

$$C_{leakage} [pF] = 4 \times d \times L$$

where:

d = smallest outer dimension of the insulation, in mm, and

L = length of sample insulation immersed in saline bath, in cm;

- *HF dielectric strength test according to 201.8.8.3.103 with an HF test voltage of 500 V_{peak} . No breakdown of the insulation shall occur;*
- *mains frequency dielectric strength test according to 201.8.8.3.104 with a test voltage of 2 100 V_{peak} . No breakdown of the insulation shall occur.*

Customer: Bradley Hyde- No. of User(s): 1 - Company: Latham and Watkins

Order No.: WS-2023-008299 - IMPORTANT: This file is copyright of IEC, Geneva, Switzerland. All rights reserved.

This file is subject to a licence agreement. Enquiries to Email: sales@iec.ch - Tel.: +41 22 919 02 11

201.15.101.5 * NE thermal performance

An NE shall not subject a PATIENT to a risk of thermal injury at the NE application site under conditions of NORMAL USE and when applied in accordance with instructions for use.

Compliance is checked by the following test.

For an NE with the PATIENT weight range marked as follows, the maximum temperature rise of any 1 cm square area under and extending 1 cm beyond the NE contact site on a PATIENT shall not exceed 6 °C immediately after a 60 s application of the specified test current, I_{test} .

Table 201.103 – Test currents by weight range

PATIENT weight range	I_{test} mA
< 5 kg	350
5 kg to 15 kg	500
> 15 kg or unspecified	700

For all MONITORING NE the contact area shall be A_a , the alarm area, as evaluated in the compliance test for subclause 201.8.4.101.

For all other NE the contact area shall be the area of the NE when applied according to the instructions for use.

For NES intended for use on small PATIENTS, these tests may be performed on adult subjects. The test surface to which the NE under test is applied shall be the skin of human subjects, or electrically and thermally equivalent surrogate media or test devices. These tests shall be repeated using a minimum of four different samples of the NE under test on each human subject or surrogate media. Where a surrogate medium or test device is used, at least 10 different samples of the NE shall be tested.

The NE and test surface temperatures of surrogate media or test devices shall be 23 ± 2 °C, and a reference temperature scan of the test surface shall be recorded immediately prior to application of the NE to the test surface. The NE shall be applied to the test surface in accordance with supplied instructions for use, except that contact area shall be A_a . The NE shall rest on the test surface for 30 min in a stable temperature environment before the application of the test current. If a thermally equivalent surrogate medium or test device is used the test may commence once thermal equilibrium is achieved.

The test current, I_{test} , applied to the electrode under test shall have an approximately sinusoidal HF waveform, and must be attained within 5 s of the beginning of the test and maintained between 100 % and 110 % of I_{test} for $60 \text{ s} \pm 1 \text{ s}$.

A second temperature scan of the test surface shall be completed within 15 s following cessation of the test current. Upon comparison with the reference scan, the temperature rise of any 1 cm square area shall not exceed 6 °C.

The temperature scanning apparatus shall have an accuracy of better than 0,5 °C and a spatial resolution of at least one sample per square cm over the entire NE contact area plus the area extending 1 cm beyond the edge of that area. Spatial correlation between the reference and second temperature scans shall be within $\pm 1,0 \text{ cm}$.

Where human subjects are employed, they shall comprise a pool of at least five males and five females having a variety of skin tissue morphologies, i.e. thin, average and thick layers of subcutaneous body fat.

Customer: Bradley Hyde- No. of User(s): 1 - Company: Latham and Watkins

Order No.: WS-2023-008299 - IMPORTANT: This file is copyright of IEC, Geneva, Switzerland. All rights reserved.

This file is subject to a licence agreement. Enquiries to Email: sales@iec.ch - Tel.: +41 22 919 02 11

Any surrogate medium or test device shall bear documented evidence that it is expected to yield temperature rise results no smaller than those from this test protocol as applied to at least 20 human subjects.

201.15.101.6 * NE contact impedance

The impedance of the electrical contact between the surface of the NE application site and the NE cord connection shall be low enough to prevent a risk of PATIENT burn due to ohmic heating during passage of HF surgical current.

For conductive NE, contact impedance shall not exceed 50 Ω , and for capacitive NES, contact capacitance shall be no less than 4 nF over the frequency range of 200 kHz to 5 MHz.

NOTE For purposes of this standard, unless otherwise specified by the manufacturer, a conductive NE presents a contact impedance with a phase angle of less than 45° at 200 kHz, and a capacitive NE a 200 kHz phase angle of 45° or greater.

Compliance is checked by the following test using at least 10 random samples of the NE under test.

The NE under test is placed in full and firm contact on a flat metallic plate. A true r.m.s. responding a.c. voltmeter having an input impedance greater than 2 k Ω and accuracy of better than 5 % over the 200 kHz – 5 MHz range is connected between the plate and the NE cord conductors in order to measure voltage U_{test} . An essentially sinusoidal test current, I_{test} , of approximately 200 mA and frequency f_{test} in the range of 200 kHz – 5 MHz is passed between the NE cord and the plate and monitored by use of a suitable true r.m.s. a.c. ammeter.

U_{test} and I_{test} are recorded at $f = 200$ kHz, 500 kHz, 1 MHz, 2 MHz and 5 MHz. For each f_{test} , contact impedance Z_c is computed as:

$$Z_c = \frac{U_{test}}{I_{test}}$$

and contact capacitance C_c is computed as:

$$C_c [nF] = \frac{I_{test} \times 10^6}{2\pi \times f_{test} \times U_{test}}$$

where

I_{test} is the r.m.s. HF test current in A;

U_{test} is the r.m.s. HF test voltage in V;

f_{test} is the HF test voltage frequency in kHz.

201.15.101.7 * NE adhesion

For NES, except MONITORING NES and NES marked for use with PATIENTS weighing less than 15 kg, if the instructions for use indicate that the NE is adhesively attached to the PATIENT, the peel strength of the adhesive shall be adequate to ensure a safe degree of contact under expected conditions of use.

Compliance is checked by the following tests.

Customer: Bradley Hyde- No. of User(s): 1 - Company: Latham and Watkins

Order No.: WS-2023-008299 - IMPORTANT: This file is copyright of IEC, Geneva, Switzerland. All rights reserved.

This file is subject to a licence agreement. Enquiries to Email: sales@iec.ch - Tel.: +41 22 919 02 11

For NES intended for use on small PATIENTS, these tests may be performed on adult subjects. Surrogate test surfaces that are shown to be equivalent to human subjects may be used.

a) Pull test

At least two samples of the NE under test are applied to convenient locations on at least 10 male and 10 female human subjects, according to instructions for use. After application, NES are allowed to remain undisturbed for 5 min to 10 min. For NES intended for use on adult PATIENTS, the attached NE cord is subjected for 10 min to a 10 N force directed along each of two orthogonal axes in a plane parallel to the skin surface at the NE cord connection point. One of the axes shall consist of the minor dimension of the NE at that point. No more than 5 % of the NE adhesive area shall separate from the skin surface in at least 90 % of the tests.

b) Conformability test

NES under test are applied to at least 5 male and 5 female human subjects on approximately cylindrical sites (e.g., extremities) having circumferences from 1,0 to 1,25 times the length of the major axis of the NE, with the major axis of the NE encircling the site. No more than 10 % of the adhesive area of the NE shall have separated from the skin surface at 1 h after application.

NOTE The conformability test is not required where this kind of application site is counter indicated in the instructions for use.

c) Fluid tolerance test

The NES are placed on at least 5 male and 5 female human subjects. The appropriate connector is connected to the NE if the NE is intended for use with a reusable cable. One litre of 0,9 % saline is poured for 5 s to 15 s from a height of 300 mm directly over the NE. No more than 10 % of the adhesive area of the NE shall have separated from the skin surface within 15 min after the saline is poured.

201.15.101.8 * NE shelf life

NES marked for single use shall comply with the requirements of 201.15.101.5 through 201.15.101.7 on the expiration date specified by the NE MANUFACTURER. Test samples may be produced by actual storage of the NES according to their instructions for use, or by accelerated aging of the NES through a cycle which has been shown to be at least as severe as equivalent recommended storage condition aging.

Compliance shall be verified by testing devices within 30 days of the expiration date or the date when accelerated aging is completed.

201.16 ME SYSTEMS

Clause 16 of the general standard applies.

201.17 Electromagnetic compatibility of ME EQUIPMENT and ME SYSTEMS

Clause 17 of the general standard applies.

202 * Electromagnetic compatibility – Requirements and tests

IEC 60601-1-2:2007 applies except as follows:

202.6.1 EMISSIONS**202.6.1.1.1 Requirements**

Addition:

aa) These requirements and tests do not apply to HF SURGICAL EQUIPMENT when the HF output is energized.

bb) HF SURGICAL EQUIPMENT shall comply with the requirements of CISPR 11 group 1, when it is switched on and in an idle state with the HF output not energized. The MANUFACTURER shall declare whether the HF SURGICAL EQUIPMENT is Class A or Class B according to its intended use.

202.6.2 IMMUNITY**202.6.2.1.10 Compliance criteria**

Addition:

The following shall be DEGRADATION of performance that does not affect ESSENTIAL PERFORMANCE or BASIC SAFETY:

- the interruption of HF power output or reset into standby mode when clearly indicated on the operation panel of HF SURGICAL EQUIPMENT;
- a change in the delivered output power as allowed in 201.12.1.101.

Compliance shall be considered to be met if the requirements of IEC 60601-1-2 are met with the above changes.

208 General requirements, tests and guidance for alarm systems in medical electrical equipment and medical electrical systems

IEC 60601-1-8:2006 applies except as follows:

Amendment:

The audible alarm and the red indicator warning light described in 201.8.4.101 shall not be considered an ALARM SIGNAL as defined in this collateral standard.

The audible signal described in 201.12.4.2.101 shall not be considered an ALARM SIGNAL as defined in this collateral standard.

Annexes

The annexes of the general standard apply.

Annex AA (informative)

Particular guidance and rationale

AA.1 General guidance

This annex provides a concise rationale for the important requirements of this particular standard and is intended for those who are familiar with the subject of the standard but who have not participated in its development. An understanding of the reasons for the main requirements is considered to be essential for the proper application of the standard. Furthermore, as clinical practice and technology change, it is believed that a rationale for the present requirements will facilitate any revision of the standard necessitated by these developments.

NOTE Testing of these devices for compliance or operation when the HF is turned on may cause test equipment to operate outside of its normal operation due to the HF electric field exposure. Suitable precautions and checks of the test instrumentation should be taken into account. This situation may also occur with the medical support instrumentation near the device.

AA.2 Rationale for particular clauses and subclauses

The following are rationales for specific clauses and subclauses in this particular standard, with clause and subclause numbers parallel to those in the body of the document.

Subclause 201.1.1 – Scope

The scope does not include equipment for cautery, i.e. for medical treatment with electrically heated metal rods or wire loops. This edition provides, to the extent feasible, separate requirements and tests for HF SURGICAL EQUIPMENT and HF SURGICAL ACCESSORIES, independent of manufacture. ASSOCIATED EQUIPMENT is included in the definition of ACCESSORIES.

Definition 201.3.207 – ASSOCIATED EQUIPMENT

Examples of ASSOCIATED EQUIPMENT are argon beam adaptors, accessory leakage monitors, NEUTRAL ELECTRODE contact monitors, and the like. The text "not intended for independent use" means the devices are used only in conjunction with HF SURGICAL EQUIPMENT. See Figure AA.1.

Definition 201.3.208 – BIPOLAR

This term is intended to apply equally to equipment and ACCESSORIES and thus is distinct from, and could possibly supplant, that of the existing subclause 201.3.209 (BIPOLAR ELECTRODE).

Definition 201.3.213 – CREST FACTOR

The measurement of CREST FACTOR is mathematically simple but difficult to carry out in a reliable manner. The r.m.s. voltage is particularly difficult to measure. The definition states that the measurements should be made in an open circuit condition. This means that the normal loads seen on the output of HF SURGICAL EQUIPMENT are not present. The load presented by the high voltage probe used to measure these voltages (10 MΩ to 100 MΩ typical) is considered to be essentially an open circuit. The following is a suggested method for these measurements that has shown a reasonable accuracy.

Customer: Bradley Hyde- No. of User(s): 1 - Company: Latham and Watkins

Order No.: WS-2023-008299 - IMPORTANT: This file is copyright of IEC, Geneva, Switzerland. All rights reserved.

This file is subject to a licence agreement. Enquiries to Email: sales@iec.ch - Tel.: +41 22 919 02 11

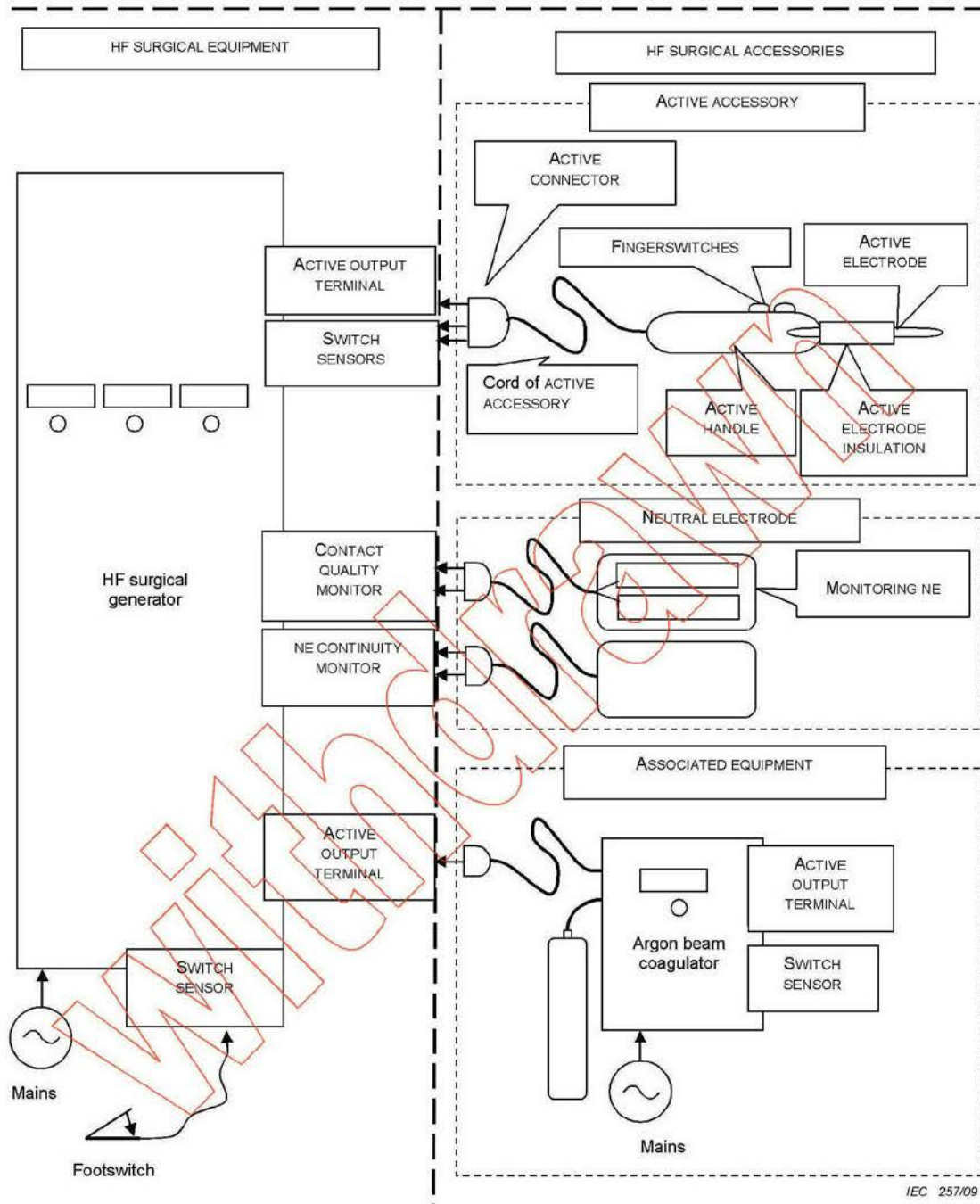


Figure AA.1 – Example of various parts of an HF surgical system

The measurements should be made from the output to the NE for MONOPOLAR outputs and across the two output poles for BIPOLAR outputs using a 1 000x or 100x high voltage probe connected to a high quality digital storage oscilloscope (DSO) with automatic measurement capabilities.

First the exact period of the signal is then measured. For continuous sinusoidal waveforms ($cf = 1,4$) that is the reciprocal of the fundamental frequency of the waveform. For non-continuous waveforms, the time period of the bursts is measured. For example, a COAGULATION waveform may have a fundamental frequency of 400 kHz with a burst repetition rate of 20 kHz. It is the precise measurement of the 20 kHz burst repetition rate that is needed. Once this time period is measured, the time base of the DSO should be modified to make the entire screen hold between 5 and 10 exact periods. For example, if the burst repetition rate is exactly 20 kHz, the period will be 50 μ s. By setting the time base of the DSO to 50 μ s per division, you should get exactly 10 waveform bursts on the screen.

The waveform is then captured and stored. Measure and record the MAXIMUM OUTPUT VOLTAGE (the absolute value of the largest peak). The r.m.s. voltage is then calculated. The most reliable method is to set the DSO to calculate the r.m.s. of the entire screen. Since the time base was adjusted to capture an exact multiple of waveforms, the calculation of r.m.s. voltage should be accurate.

An alternate way of measuring the r.m.s. voltage would be to connect the output of the high voltage probe into a thermal sensing true r.m.s. voltmeter that is RATED for the CREST FACTOR of the waveform being measured.

The CREST FACTOR can now be calculated.

Definition 201.3.214 – CUTTING

It is generally believed that HF surgical CUTTING involves microscopic cellular ablation resulting from short electrical sparks being struck between the ACTIVE ELECTRODE and the tissue.

Definition 201.3.215 – EARTH REFERENCED PATIENT CIRCUIT

For purposes of this particular standard, the impedance of this path, at the lowest HF operating frequency, is 10 Ω or less. See Figure AA.3.

Definition 201.3.217 – FULGURATION

FULGURATION generally requires HF peak output voltages of at least 2 kV in order to ignite and sustain the long sparks. This mode is also known as spray or non-contact COAGULATION and may be enhanced by incorporation of a stream of inert gas such as argon.

Definition 201.3.218 – HIGH FREQUENCY

Frequencies above 200 kHz should be used for MONOPOLAR applications in order to avoid the unwanted stimulation of nerves and muscles which would result from the use of low frequency current. Lower frequencies may be used for BIPOLAR techniques if the RISK ANALYSIS shows the possibility of neuromuscular stimulation has been mitigated to an acceptable level.

Normally, frequencies above 5 MHz are not used in order to minimize the problems associated with HIGH FREQUENCY LEAKAGE CURRENTS. However, higher frequencies may be used in the case of BIPOLAR techniques. It is generally recognized that 10 mA is the lower threshold of thermal effects on tissue.

Customer: Bradley Hyde- No. of User(s): 1 - Company: Latham and Watkins

Order No.: WS-2023-008299 - IMPORTANT: This file is copyright of IEC, Geneva, Switzerland. All rights reserved.

This file is subject to a licence agreement. Enquiries to Email: sales@iec.ch - Tel.: +41 22 919 02 11

Definition 201.3.223 – HF SURGICAL MODE

The term HF SURGICAL MODE should be clearly distinct from “mode of operation” as used in subclauses 6.6 and 7.2.11 of the general standard in reference to operational duty cycle.

Definition 201.3.224 – MAXIMUM OUTPUT VOLTAGE

This parameter is intended for comparison by the OPERATOR to RATED ACCESSORY VOLTAGE to ensure safety.

Definition 201.3.225 – MONITORING NE

A CONTACT QUALITY MONITOR is functional only when used with a MONITORING NE. A MONITORING NE is also known as a split or divided plate since the conductive area is split into two or more parts.

Definition 201.3.226 – MONOPOLAR

This definition is intended to apply equally to equipment and ACCESSORIES and thus is distinct from that of the existing subclause 201.3.203 (ACTIVE ELECTRODE).

Subclause 201.4.2 – RISK MANAGEMENT PROCESS for ME EQUIPMENT or ME SYSTEMS

In MONOPOLAR surgery, the three elements that are used as a system are the HF SURGICAL ACCESSORY, the HF generator and the NEUTRAL ELECTRODE. MANUFACTURERS of any one or more of these elements need to consider the possible use of their products during non-traditional, high current situations. These situations might include, but are not limited to: tissue lesioning, tissue ablation, tissue vaporization, and procedures where conductive fluid is introduced into the surgical site for distension or to conduct the HF current. In high current situations, there is a risk that heating under a NEUTRAL ELECTRODE may be high enough to cause HARM to the PATIENT.

Subclause 201.4.3 – ESSENTIAL PERFORMANCE

After a careful consideration of the clauses within this standard, it was decided that they all deal with BASIC SAFETY as defined in the general standard. Manufacturers have the ability to identify functions of HF SURGICAL EQUIPMENT which are considered ESSENTIAL PERFORMANCE in accordance with their RISK MANAGEMENT PROCESS.

Subclause 201.5.4 – Other conditions

Instruments used to measure HF currents, including HF voltmeter/current sensor combinations, should register true r.m.s. with a total verified accuracy of 5 % of reading or better from 10 kHz to at least 5 times the fundamental frequency of the HF SURGICAL MODE being tested. HF output instruments should register to specified accuracy within 3 s of application of the measured variable. Transient readings of less than 1 s duration may be ignored.

Resistors used for HF testing should be RATED at no less than 50 % of the expected power dissipation expected for a given test, and should present a resistive component of impedance within 3 % of the specified value and no more than 8,5 degrees of impedance phase from 10 kHz to 5 times the fundamental frequency for the HF SURGICAL MODE being tested.

Instruments used for measuring HF voltages should be RATED at no less than 150 % of the expected peak voltage and should have verified accuracy of 5 % of reading or better from 10 kHz to 5 times the fundamental frequency of the signal being measured.

For each HF SURGICAL MODE, the term “fundamental frequency” means the frequency of the highest amplitude spectral line of the measured HF output voltage when operated at maximum power setting into an open circuit.

Customer: Bradley Hyde- No. of User(s): 1 - Company: Latham and Watkins

Order No.: WS-2023-008299 - IMPORTANT: This file is copyright of IEC, Geneva, Switzerland. All rights reserved.

This file is subject to a licence agreement. Enquiries to Email: sales@iec.ch - Tel.: +41 22 919 02 11

A primary objective of the revision of this particular standard is to isolate HF ACCESSORY requirements and tests from any specific HF SURGICAL EQUIPMENT. Further, this standard should clearly specify instrumentation for required tests to ensure repeatability of results, particularly for test agencies which may not be conversant with accepted HF test methods. Due to the brevity of power application and the greater availability of lower-power resistors which satisfy the low reactance requirement, resistors RATED as low as 50 % of expected power are suitable, but no lower.

Subclause 201.7.4.2 – Control devices

As the power delivered to the load depends on the load resistance, a graduation in relative units is considered to be adequate. However, if an output indication displays the actual power output in watts, it must do so over the total range of load resistance, otherwise the power delivered to the PATIENT may differ from that indicated and hence create an unacceptable RISK. If the numeral "0" is displayed, the OPERATOR will expect zero output at this position of the control.

Subclause 201.7.8.1 – Colours of indicator lights

The standardization of the colours of indicator lights is regarded as a safety feature.

For many years the yellow indicator light has been used to signify that the CUTTING mode is selected or in use on HF SURGICAL EQUIPMENT. During surgery, a "blend" mode is used mainly for CUTTING with varying amounts of COAGULATION added. As the main function of "blend" is to cut, it is considered that a yellow light is most appropriate when "blend" is in use.

Subclause 201.7.8.2 – Colours of controls

The same colour coding as specified for indicator lights should be used in other places to avoid confusion.

Subclause 201.7.9.2.2.101 a)

The advice concerning avoidance of unwanted burns is based on experience. In particular:

- 1) In past editions of this standard this section included advice to place the NEUTRAL ELECTRODE as close to the operating field as possible. In general, minimising the distance between the operating field and the NEUTRAL ELECTRODE reduces the load resistance and, for a given power at the site of the ACTIVE ELECTRODE, the power output required from HF SURGICAL EQUIPMENT and also the HF voltage across the PATIENT. However, if the direct path between the ACTIVE ELECTRODE and the NEUTRAL ELECTRODE includes small cross sectional areas of tissue, the current density could cause undesired heating and tissue damage. Therefore the OPERATOR should rely on the instructions for use provided by the MANUFACTURER of the NEUTRAL ELECTRODE for specific placement instructions.
- 2) Small area contacts with objects having a low impedance to earth at HIGH FREQUENCIES may result in high current densities and hence unwanted burns.
- 3) There may be some HF voltage difference between these parts of the PATIENT's body which may cause an unwanted current to flow.
- 4) The current flowing to the leads of the monitoring equipment may cause burns at the site of the monitoring electrodes.
- 5) The capacitance between the electrode cable and the PATIENT may result in some local high current densities.
- 6) In certain cases, BIPOLAR technique can avoid unwanted tissue damage, especially where bony structures having a relatively high resistance or parts of the body having a relatively small cross section are involved.
- 8) In this case, the application of the NEUTRAL ELECTRODE and its connections should be checked before selecting a higher output power.

Customer: Bradley Hyde- No. of User(s): 1 - Company: Latham and Watkins

Order No.: WS-2023-008299 - IMPORTANT: This file is copyright of IEC, Geneva, Switzerland. All rights reserved.

This file is subject to a licence agreement. Enquiries to Email: sales@iec.ch - Tel.: +41 22 919 02 11

Not all advice is necessary, if only a BIPOLAR output or a RATED OUTPUT POWER not exceeding 50 W without NEUTRAL ELECTRODE is available.

Subclause 201.7.9.2.2.101 c)

IEC 60601-2-18:1996 contains requirements prescribing that MANUFACTURERS of HF energized endotherapy devices shall specify them as suitable for a certain maximum allowed HF peak output voltage which shall be defined in the ACCOMPANYING DOCUMENTS for such ACCESSORIES.

However IEC 60601-2-18-A1 (2000) requires a "RATED recurring peak voltage" together with the "mode(s) of intended use". The experts of Maintenance Team 17 felt that this information on one hand is insufficient, as the modes of intended use like "spray COAGULATION" are not clearly technically defined and may vary much between different brands and models of HF SURGICAL EQUIPMENT. On the other hand, it was considered as impracticable to give such rather complex information to the user of the equipment.

Therefore it was considered as more practicable to provide the user only with a RATED ACCESSORY VOLTAGE and a MAXIMUM OUTPUT VOLTAGE for any output setting in order to enable the user to judge whether any HF SURGICAL ACCESSORY or ASSOCIATED EQUIPMENT can be safely used with any certain output setting of the generator.

At HIGH FREQUENCY the stability of insulation is affected by dielectric heating so the relationship between the MAXIMUM OUTPUT VOLTAGE and the CREST FACTOR is important.

Further, it was considered that with all currently known brands and models of generators in modes and settings producing higher output voltages, the CREST FACTOR is always increased along with the voltage. Therefore a general relation between output voltage and CREST FACTOR was developed as shown in Figure AA.2.

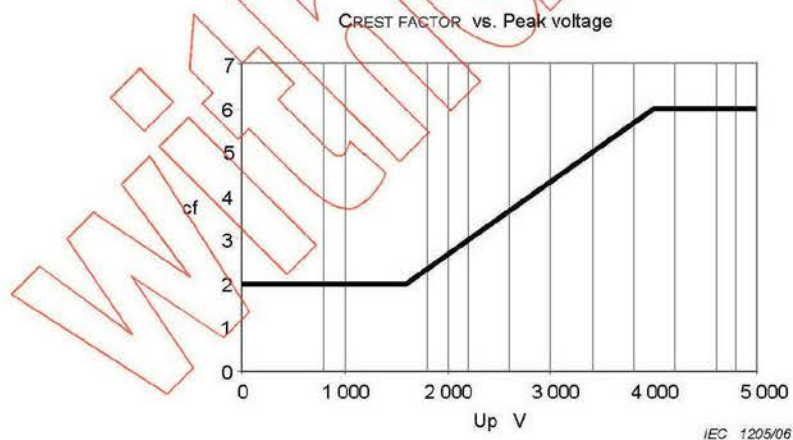


Figure AA.2 – CREST FACTOR vs. peak voltage

A safe situation exists whenever a RATED ACCESSORY VOLTAGE is matched to an output voltage of an HF SURGICAL EQUIPMENT having a CREST FACTOR which falls on or above the line in the diagram. The RATED ACCESSORY VOLTAGE must not be less than the MAXIMUM OUTPUT VOLTAGE, since the HF SURGICAL ACCESSORY or ASSOCIATED EQUIPMENT must fulfil the requirements of 201.8.8.3.103 which takes into account the CREST FACTOR.

Provision is made for the case in which a generator at a certain setting has a MAXIMUM OUTPUT VOLTAGE with a corresponding CREST FACTOR which falls below the line. In this case, to ensure safety, the RATED ACCESSORY VOLTAGE must be high enough to ensure that there is no insulation breakdown of the HF SURGICAL ACCESSORY or ASSOCIATED EQUIPMENT when used with

Customer: Bradley Hyde- No. of User(s): 1 - Company: Latham and Watkins

Order No.: WS-2023-008299 - IMPORTANT: This file is copyright of IEC, Geneva, Switzerland. All rights reserved.

This file is subject to a licence agreement. Enquiries to Email: sales@iec.ch - Tel.: +41 22 919 02 11

that particular HF SURGICAL EQUIPMENT in that particular HF SURGICAL MODE at that particular output setting. This precaution is necessary in order to take into account the dielectric heating produced by lower CREST FACTOR waveforms. The safe value of RATED ACCESSORY VOLTAGE must be found by testing the HF SURGICAL ACCESSORY or ASSOCIATED EQUIPMENT with the HF SURGICAL EQUIPMENT.

Subclause 201.7.9.2.2.101 e)

The OPERATOR must know which MONITORING NES are safe and functional with the CQM. Many OPERATORS mistakenly believe that with the advent of CQM, intraoperative surveillance of NE contact is no longer necessary.

Subclause 201.7.9.2.2.101 g)

Although the required measures in 201.8.4.102 are intended to reduce neuromuscular stimulation significantly, it cannot be completely eliminated especially when electrical arcs are produced. Therefore a warning is necessary to make the user aware that in sensitive structures neuromuscular stimulation can still occur leading to secondary risks like injury caused by muscle contractions. See also the rationale for 201.8.4.102.

Subclause 201.7.9.2.2.101 h)

For systems of HF surgical devices used under these conditions, there is an increased concern for NEUTRAL ELECTRODE burns.

Subclause 201.7.9.2.14 – Accessories, supplementary equipment, used material

Some OPERATORS believe incorrectly that CQM is intrinsic to either the CONTACT QUALITY MONITOR or MONITORING NE alone. It is important that all OPERATORS understand all of the physical requirements necessary to achieve CQM functionality.

Subclause 201.7.9.2.14 e)

This information should enable the OPERATOR to judge the suitability of an HF SURGICAL EQUIPMENT or its output setting for a particular ACCESSORY with regard to its isolation quality.

Subclause 201.7.9.2.14 f)

The OPERATOR must know which CQM(s) are operative with a given NE.

Subclause 201.7.9.2.14 g)

A statement of compatibility may take different forms as long as it can be understood by the OPERATOR (e.g. an impedance based CQM system where the alarm sounds based on the following conditions ..., a CQM system found in the following list of equipment ..., a CQM system from the following MANUFACTURERS ..., as well as other forms).

Subclause 201.7.9.3.1 – General

Some specialized HF SURGICAL EQUIPMENT do not provide OPERATOR adjustable output settings.

These diagrams should enable the OPERATOR to judge the suitability of an HF SURGICAL EQUIPMENT for a particular purpose. If the HF SURGICAL EQUIPMENT has discreet blend selections (e.g. blend 1, blend 2, etc), then a diagram would be created for each discreet mode. If the HF SURGICAL EQUIPMENT has a variable blend control where the setting may be

Customer: Bradley Hyde- No. of User(s): 1 - Company: Latham and Watkins

Order No.: WS-2023-008299 - IMPORTANT: This file is copyright of IEC, Geneva, Switzerland. All rights reserved.

This file is subject to a licence agreement. Enquiries to Email: sales@iec.ch - Tel.: +41 22 919 02 11

continuously adjusted, then the control should be set to the blend setting that provides the greatest hemostasis.

Subclause 201.8.4.101 – NEUTRAL ELECTRODE monitoring circuit

Undetected interruption of the NEUTRAL ELECTRODE cable in HF SURGICAL EQUIPMENT or insufficient electrical contact between the NEUTRAL ELECTRODE and the PATIENT may lead to severe burns. Therefore, as a minimum requirement, monitoring of a failure of the NEUTRAL ELECTRODE circuit or its connections is required for such HF SURGICAL EQUIPMENT having a RATED OUTPUT POWER in excess of 50 W.

The revised subclause title is intended to distinguish over the various other monitor circuits which may be present in HF SURGICAL EQUIPMENT, such as output power fault detection, and the like.

A CONTACT QUALITY MONITOR should be shown to function effectively when used with any MONITORING NE listed as compatible. When combined with new requirements for NE thermal performance, the risk of NE site burns is effectively mitigated. Because of the technical variations and proprietary nature of existing CQM schemes, imposition of a fully ACCESSORY independent requirement is judged impractical.

Full contact means that the NE has been applied according to the instructions for use such that the conductive portion within the NE is as close to the human subject (or suitable surrogate surface) as possible without any voids or spaces.

The literature listed as items [1] through [4]³⁾ in the bibliography is recommended as a guide in evaluating suitable surrogate surfaces.

Subclause 201.8.4.102 – Neuromuscular stimulation

Due to the rectifying effect of arcs between the ACTIVE ELECTRODE and tissue, d.c. and low frequency components may cause neuromuscular stimulation. This undesirable stimulation is effectively reduced by the use of an appropriate value of series capacitance and shunt resistance.

Subclause 201.8.5.1.2 – MEANS OF PATIENT PROTECTION (MOPP)

These reduced requirements are considered to be adequate because the "voltages stressing the insulation..." are of HIGH FREQUENCY and therefore, if insulation fails between the HF APPLIED PARTS and the ENCLOSURE, the RISK is much lower than at lower frequencies.

Subclause 201.8.5.2.3 – PATIENT leads

This subclause of the general standard is designed to prevent a connection between the PATIENT and either ground or a hazardous voltage. This subclause assumes that a connection may occur at any time, and that the contact with the PATIENT is either continuous or unsupervised.

The situation with HF SURGICAL ACCESSORIES is quite different, because this kind of equipment is intended to be used only under the control of a doctor or trained medical staff. Possible hazardous situations, which may occur by insertion of connectors of NEUTRAL ELECTRODES into mains connectors, such as mains outlets or sockets of detachable power supply cords, are covered by this subclause of the particular standard.

³⁾ Figures in square brackets refer to the Bibliography.

Unlike electrocardiographic monitoring electrodes which can be expected to be applied by OPERATORS untrained in electrical HAZARDS, HF SURGICAL EQUIPMENT and ACCESSORIES are accessible only to OPERATORS highly qualified and trained in restricted access locations.

ACTIVE and BIPOLAR ELECTRODES are applied only under the direct control of a surgeon who may be expected to interrupt contact with the PATIENT at the slightest sign of an unexpected PATIENT's response.

Subclause 201.8.6.1 – Applicability of requirements

For low powered MONOPOLAR HF SURGICAL EQUIPMENT used without a NEUTRAL ELECTRODE this is common practice; it is considered not to create any safety problem.

Subclause 201.8.7.1 – General requirements

The requirements for LEAKAGE CURRENT specified in the general standard are intended to provide protection against the risk of electric shock.

In this particular standard some requirements for HF LEAKAGE CURRENT are also given in order to reduce the risk of unwanted burns.

This subclause of the general standard is concerned with LEAKAGE CURRENTS which lead to electric shock, not with therapeutic currents such as are produced by HF SURGICAL EQUIPMENT. Appropriate tests for HF LEAKAGE CURRENT for HF SURGICAL EQUIPMENT with multiple PATIENT circuits are given in 201.8.7.3.101 c) cross-coupling between different HF PATIENT CIRCUITS.

Subclause 201.8.7.3 – Allowable values


Monitoring currents which flow exclusively between the parts of a split NEUTRAL ELECTRODE are considered not to need limitation according to TYPE CF APPLIED PARTS, independent of the degree of protection against electric shock (TYPE BF or CF APPLIED PARTS), because these currents can be expected never to flow across the heart.

Subclause 201.8.7.3.101 a) – Thermal effects of HF LEAKAGE CURRENTS

HF SURGICAL EQUIPMENT designed for use without a NEUTRAL ELECTRODE had to be exempted since, in such HF SURGICAL EQUIPMENT, a differentiation between functional and HF LEAKAGE CURRENT is impossible. Therefore, the measurement of functional and HF LEAKAGE CURRENT is meaningless.

As distinct from the LEAKAGE CURRENT measurements of the general standard, a measuring resistance of 200 Ω is specified here to simulate the load impedances prevailing in actual situations so as to give the maximum leakage power. The values specified result in a power of 2 W, which is considered to be a reasonable limit. Test 2 of the earth referenced case is specified to verify that the impedance to earth at HIGH FREQUENCY is sufficiently low.

An earthed conductive plane under the insulating table and bundling the power supply cord rather than coiling it, improve the reproducibility of the measurement considerably.

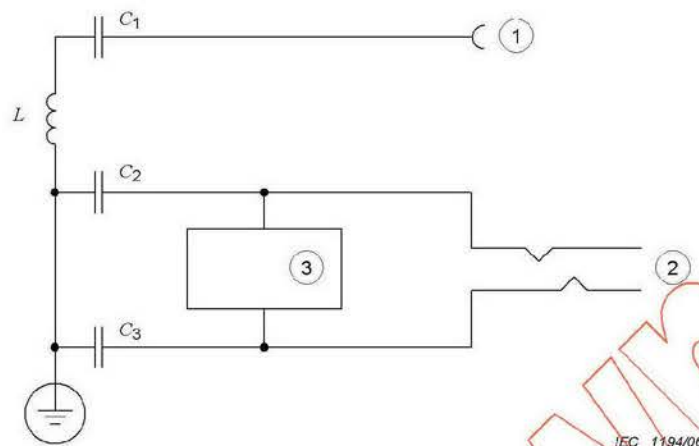**Key**

- ① Connector for ACTIVE ELECTRODE
- ② Connector for NEUTRAL ELECTRODE
- ③ Monitor

C_1 not to exceed 5 nF

$C_2 = C_3$ not to exceed 25 nF

X_{C2} and X_{C3} at operating frequency each not to exceed 20 Ω

Z_L at 50 Hz not to exceed 1 Ω

Figure AA.3 – Example of PATIENT circuit with NEUTRAL ELECTRODE referenced to earth at operating frequencies

Subclause 201.8.7.3.101 a) 3)

Experience in testing BIPOLAR HF SURGICAL EQUIPMENT has shown that these limits are reasonable and the test realistic. The RISK MANAGEMENT FILE may be reviewed for adequate explanation of alternate means of measurement and/or risk mitigation.

Subclause 201.8.7.3.101 b) – HIGH FREQUENCY LEAKAGE CURRENTS measured directly at the HF SURGICAL EQUIPMENT terminals

A test of the isolation of HF SURGICAL EQUIPMENT at HIGH FREQUENCY is easily achieved by placing load resistances and measuring devices directly on the output terminals. In this case a limit of 100 mA is specified because the contribution from the leads is not included. However, in order to ensure that all complex impedances resulting from leads and ACCESSORIES (for example ACTIVE ELECTRODES with FINGERSWITCHES) are considered, the test in 201.8.7.3.101 a) is also included.

Subclause 201.8.8.3.101 – ACTIVE ACCESSORY insulation

HF SURGICAL EQUIPMENT is capable of producing high voltages which will appear on the insulated conductive parts of HF SURGICAL ACCESSORIES. The insulation on these ACCESSORIES must withstand this voltage stress and limit the HF LEAKAGE CURRENT density appearing on exposed surfaces in order to mitigate the risk of unintended burns to the PATIENT and the OPERATOR. This insulation is subjected to considerable stress in practical use, and therefore the requirements contain a safety margin. The insulation applied to any part of an ACTIVE ACCESSORY must maintain adequate dielectric strength after extended exposure to conductive fluids and, except for accessories intended for single use, repeated sterilization.

Customer: Bradley Hyde- No. of User(s): 1 - Company: Latham and Watkins

Order No.: WS-2023-008299 - IMPORTANT: This file is copyright of IEC, Geneva, Switzerland. All rights reserved.

This file is subject to a licence agreement. Enquiries to Email: sales@iec.ch - Tel.: +41 22 919 02 11

NOTE This subclause has been completely redrafted to cover only dielectric strength of the various parts of ACTIVE ACCESSORY insulation, independent from any particular HF SURGICAL EQUIPMENT. The revised requirements and compliance tests have drawn upon the current editions of ANSI/AAMI HF18 and IEC 60601-2-18 with a goal of harmonization.

Requirements for NES are now compiled under 201.15.101.

Subclause 201.8.8.3.102 – ACTIVE ACCESSORY HF leakage

The HF leakage requirements are based on ANSI/AAMI HF18:2001, subclause 4.2.5.2. The rationale for these requirements is excerpted below. In order to use common SI units, the text and formulas for both the normative language and the rationale have been changed from the original.

The 1 MHz maximum operating frequency and the RATED ACCESSORY VOLTAGE constitute a reasonable margin between the test limits and the performance of present-day cables while maintaining a considerable margin between the test limits and that which would produce current densities of 100 mA/cm².

All of the selected values in combination permit an equivalent current density of 11,46 mA/cm², which is nearly an order of magnitude below the recognized burn threshold of 100 mA/cm² for 10 s. Therefore, while it may be argued that the levels of one or more of the factors may be higher under extreme clinical conditions, the safety margin built into the requirements is judged to be sufficient.

The cables of NEUTRAL ELECTRODES are allowed twice the leakage of the cables of ACTIVE ACCESSORIES because the voltage levels developed between the conductors of such cables and the PATIENT's skin are generally much lower. BIPOLAR accessories are allowed twice the leakage of MONOPOLAR cables because the voltage of use is generally much lower than with the MONOPOLAR mode.

The following allowances have been incorporated in this particular standard to permit use of ordinary HF SURGICAL EQUIPMENT to generate test voltages:

The allowable test voltage range for MONOPOLAR ACCESSORIES should exceed the Paschen minimum of about 280 V_{peak} in order to permit corona development, but need not exceed the typical CUT output voltage of about 1 000 V_{peak}. Nor should the peak test voltage exceed the RATED ACCESSORY VOLTAGE.

These allowances are accommodated in harmonization with ANSI/AAMI HF18 by adjusting the HF LEAKAGE CURRENT conformance limit as follows:

$$I_{\text{leakage}} [\text{mA}] = 9,0 \times 10^{-6} \times d \times L \times f_{\text{test}} \times U_{\text{peak}}$$

For BIPOLAR cords and NE cords, the HF LEAKAGE CURRENT is doubled:

$$I_{\text{leakage}} [\text{mA}] = 1,8 \times 10^{-5} \times d \times L \times f_{\text{test}} \times U_{\text{peak}}$$

The risk of HF LEAKAGE CURRENT passing through ACTIVE ELECTRODE INSULATION and the insulation of NE cords is deemed at least as serious as that of cords of ACTIVE ACCESSORIES, and therefore those parts are included in these requirements.

Alternate HF LEAKAGE test:

The equivalent capacitance of the ANSI/AAMI HF18 HF LEAKAGE test path is derived as follows:

Customer: Bradley Hyde- No. of User(s): 1 - Company: Latham and Watkins

Order No.: WS-2023-008299 - IMPORTANT: This file is copyright of IEC, Geneva, Switzerland. All rights reserved.

This file is subject to a licence agreement. Enquiries to Email: sales@iec.ch - Tel.: +41 22 919 02 11

Given

$$I_{\text{leakage}} [\text{A}] = \frac{U_{\text{test}} [\text{V}]}{X_{\text{leakage}} [\Omega]}$$

and

$$X_{\text{leakage}} [\Omega] = \frac{1}{(2\pi \times f_{\text{test}} [\text{Hz}] \times C [\text{F}])}$$

then

$$I_{\text{leakage}} [\text{mA}] \times 10^{-3} = U_{\text{test}} [\text{V}] \times f_{\text{test}} [\text{kHz}] \times 10^3 \times 2\pi \times C [\text{pF}] \times 10^{-12}$$

thus

$$C [\text{pF}] = \frac{I_{\text{leakage}} [\text{mA}] \times 10^6}{[2\pi \times U_{\text{test}} [\text{V}] \times f_{\text{test}} [\text{kHz}]]} \quad (\text{AA.1})$$

The r.m.s. value of a sinusoidal test voltage is evaluated as:

$$V = \frac{V_{\text{p-p}}}{2\sqrt{2}} = 0,3536 \times V_{\text{p-p}}$$

The constants used for the HF leakage test are:

$$V_{\text{p-p}} = 800 [\text{V}];$$

$$U_{\text{test}} = 282,8 [\text{V}];$$

$$f_{\text{test}} = 1\,000 [\text{kHz}];$$

$$I_{\text{leakage}} = 3,6d \times L [\text{mA}].$$

The limiting capacitance according to Equation AA.1 is thus:

$$C [\text{pF}] = 2,026 \times d [\text{mm}] \times L [\text{cm}]$$

for all but BIPOLAR ACTIVE ACCESSORIES and NE cords. These are allowed twice the LEAKAGE CURRENT which yields:

$$C [\text{pF}] = 4,052 \times d [\text{mm}] \times L [\text{cm}].$$

For purposes of this standard, these results are rounded down to $2 \times d \times L$ and $4 \times d \times L$ [in pF] respectively.

The technical equivalence of the foregoing alternate capacitance-based test method to the precedent HF LEAKAGE CURRENT method has been validated by Keller [5] and König [6].

Subclause 201.8.8.3.103 – ACTIVE ACCESSORY HF dielectric strength

As the dielectric stress is at HIGH FREQUENCY in practice, additional testing at HIGH FREQUENCY is required. A saline test electrode reasonably simulates the wet PATIENT and OPERATOR tissue in or near the surgical site. The use of a thin wire wrapped over insulation has been shown to

Customer: Bradley Hyde- No. of User(s): 1 - Company: Latham and Watkins

Order No.: WS-2023-008299 - IMPORTANT: This file is copyright of IEC, Geneva, Switzerland. All rights reserved.

This file is subject to a licence agreement. Enquiries to Email: sales@iec.ch - Tel.: +41 22 919 02 11

induce corona discharge damage which can be detected by the subsequent mains frequency dielectric strength test. Each test was independently selected to exert worst case stress on the insulation being challenged. The measurement of V_{peak} and the CREST FACTOR should occur simultaneously with the test of the ACCESSORY to ensure that their values do not change due to loading by the accessory. During these tests, measuring the CREST FACTOR in a loaded state is acceptable.

These requirements and tests harmonize to the extent possible with IEC 60601-2-18.

Subclause 201.8.8.3.104 – ACTIVE ACCESSORY mains frequency dielectric strength

It is known that HF test voltages greater than 120 % of that available from HF SURGICAL EQUIPMENT are difficult to achieve. Step-up transformers tend to distort the HF waveform, and the capacitance of the dielectric being tested can load the HF test voltage source. In order to stress insulation with an acceptably high margin, a d.c. or mains frequency test is required. This test follows the HF dielectric strength test in order to detect any corona-induced weaknesses.

Elevated temperatures produced by dielectric stress can alter the internal structure of HF ACTIVE ACCESSORIES. Any incorporated FINGERSWITCH should function reliably and not activate its output inadvertently following all of the dielectric strength tests.

NOTE The metal foil used in the compliance test should be highly conductive.

Subclause 201.8.10.4.2 – Connection cords

The requirements of these two subclauses (derived from IEC 60601-2-4) are specified because ACTIVE ACCESSORIES and their cables are subject to considerable stress in use and typical failure modes can present a hazard to staff and/or PATIENTS. Once a cable fatigues in use, it is common that it will overheat and either ignite itself or ignite nearby materials, endangering staff and PATIENTS. These requirements will establish a reference level for durability of such cables.

Subclause 201.8.10.4.101 – SWITCH SENSORS

The output switch is required to be of a momentary type in order to prevent unintentional energization of the output. The requirement for isolated extra-low voltage takes into account the severe environmental conditions under which these footswitches, FINGERSWITCHES and their cables are used. The requirement against the effects of entry of liquids is already defined in subclause 201.11.6.5 of this particular standard.

It is considered that using one FINGERSWITCH for selecting a multiple function, for example CUTTING or COAGULATION, could result in confusion and a potential hazard if a surgeon unfamiliar with the system were to use it. One example of this is light pressure on the switch may give COAGULATION, heavier pressure may give CUTTING.

This subclause assumes the equipment is turned on.

Subclause 201.11.1.1 – Maximum temperature during NORMAL USE

The operating conditions specified here are deemed to be the most severe conditions likely to occur in practical use.

Subclause 201.11.6.3 – Spillage on ME EQUIPMENT and ME SYSTEMS

The test quantity of one litre represents a liquid filled bottle (for example an infusion solution), the presence of which in an operating room is considered to be likely.

Customer: Bradley Hyde- No. of User(s): 1 - Company: Latham and Watkins

Order No.: WS-2023-008299 - IMPORTANT: This file is copyright of IEC, Geneva, Switzerland. All rights reserved.

This file is subject to a licence agreement. Enquiries to Email: sales@iec.ch - Tel.: +41 22 919 02 11

Subclause 201.11.6.5 a) – Ingress of water or particulate matter into ME EQUIPMENT and ME SYSTEMS

A footswitch may be exposed to a considerable amount of water or other liquids during certain operations, and also when it is cleaned (for example by total immersion); consequently watertightness is required.

Revisions of the immersion test to replace inspection with functional and dielectric strength testing are under consideration. No current IEC 60529 tests are deemed appropriate for the expected operating room environment.

Subclause 201.11.6.5 b)

A certain degree of water protection has to be required for FINGERSWITCHES to prevent inadvertent activation of an output by the ingress of conductive fluids. This test is independent of specific HF SURGICAL EQUIPMENT. An a.c. impedance measurement of 1 kHz avoids measurement errors due to polarization effects in saline which may bridge the switching contacts, and the voltage is consistent with subclause 201.8.10.4.101. The impedance limit was chosen as twice the maximum threshold stipulated by 201.8.10.4.101.

Subclause 201.11.6.7 – Sterilization of ME EQUIPMENT and ME SYSTEMS

Applicable to all ACCESSORY specific requirements. The specified parts are expected to enter the sterile surgical field during use and thus will be re-sterilised after each use. There are no requirements or tests which can rationally be excepted from this requirement.

ACTIVE ACCESSORIES marked for single use are unsuitable for re-sterilisation and thus are exempted from this requirement.

Subclause 201.12.1.102 a)

In the load resistance range normally prevailing in practical use, lowering the output setting should never result in an increase in output power.

Subclause 201.12.1.102 b) Monotonicity of output control setting – BIPOLAR outputs

The RISK MANAGEMENT FILE may be reviewed for adequate explanation of alternate means of measurement.

Subclause 201.12.1.103 – Accuracy of MAXIMUM OUTPUT VOLTAGE

The maximum peak output voltage may appear at output settings other than maximum and with applied loads other than open circuit.

Subclause 201.12.2 b)

The standardization of the position of activating controls is required to reduce human errors. Controls for functions other than CUTTING and COAGULATION activation may also appear on the ACTIVE HANDLE.

Subclause 201.12.2 d)

Within this subclause, the term simultaneous activation refers to either situation described in 201.12.2 c).

Customer: Bradley Hyde- No. of User(s): 1 - Company: Latham and Watkins

Order No.: WS-2023-008299 - IMPORTANT: This file is copyright of IEC, Geneva, Switzerland. All rights reserved.

This file is subject to a licence agreement. Enquiries to Email: sales@iec.ch - Tel.: +41 22 919 02 11

In clinical use, the problems of co-ordination of the simultaneous use of more than one ACTIVE OUTPUT TERMINAL are considered to create unacceptable HAZARDS if only one output switch and set of controls are incorporated.

Subclause 201.12.2 e) and f)

Specifically places the burden for avoidance of incorrect connections on the equipment. Flying leads do not adequately discourage improper connections. Misconnection of a single pin ACCESSORY presents no conceivable HAZARD.

Examples of BIPOLAR ACTIVE CONNECTOR configurations which conform to this requirement are under consideration.

Subclause 201.12.2 g)

The pre-indication of the output and/or function (for example CUTTING or COAGULATION) is an essential safety feature where they are energized by the same output switch.

Subclause 201.12.4.2 – Indication of parameters relevant to safety

Within this subclause, the term simultaneous activation refers to the situation described in 201.12.2 c) 1).

Subclause 201.12.4.3.101 – Output reduction means

In the load resistance range normally prevailing in practical use, lowering the output setting should never result in an increase in output power.

Subclause 201.12.4.4.101 – Maximum allowed output power in SINGLE FAULT CONDITIONS

Although not required for MONOPOLAR HF SURGICAL EQUIPMENT having a RATED OUTPUT POWER not exceeding 50 W, compliance with this subclause is recommended. This requirement is intended to apply to all BIPOLAR outputs of HF SURGICAL EQUIPMENT.

Subclause 201.12.4.4.102 – Output power during simultaneous activation

Independent outputs must deliver their intended output power to prevent a hazard. This is especially true when one output is set at a level substantially lower than another, but both can be activated simultaneously.

Where multiple outputs share the power of a single mode (e.g. simultaneous COAGULATION), a hazard could exist if a single output delivers more power than the intended power or if the sum total of the power delivered in all of the simultaneously activated outputs exceeds the intended power.

Subclause 201.13.2.13.101 – Protection against the effects of short-circuiting of the electrodes

Some ACCESSORIES, for example resectoscopes or BIPOLAR ELECTRODES, may short-circuit the output in NORMAL USE and the output circuit is frequently energized while open circuited. It is considered practical to design HF SURGICAL EQUIPMENT which will not be damaged by repeated short circuiting and by the open circuiting of the output for short periods of time. The revised text is intended to eliminate a question of which BIPOLAR output terminal is the NEUTRAL ELECTRODE and whether this requirement applies to BIPOLAR outputs.

Customer: Bradley Hyde- No. of User(s): 1 - Company: Latham and Watkins

Order No.: WS-2023-008299 - IMPORTANT: This file is copyright of IEC, Geneva, Switzerland. All rights reserved.

This file is subject to a licence agreement. Enquiries to Email: sales@iec.ch - Tel.: +41 22 919 02 11

Subclause 201.15.4.1.101 and 201.15.4.1.102

The requirements of these subclauses relate to the compatibility of detachable parts of ACTIVE ACCESSORIES. This issue becomes important for third party accessories and can cause operational difficulties in clinical practice, leading to delayed or interrupted procedures.

Many ACTIVE HANDLES provide for the use of any of a variety of specialized, OPERATOR selected, detachable ACTIVE ELECTRODES. There is no standardization of electrode interface amongst the ACTIVE HANDLES of different manufacture. It is known that, although it may appear to the OPERATOR that an ACTIVE ELECTRODE from one manufacturer may fit the ACTIVE HANDLE from another, PATIENT injuries have resulted from incompatibilities such as:

- inadequate SEPARATION between the conductive parts of the ACTIVE HANDLE - ACTIVE ELECTRODE interface and PATIENT tissue;
- arcing across a gap between the intended electrical mating parts, resulting in melting and/or ignition of insulation;
- inadequate mechanical retention force, resulting in the ACTIVE ELECTRODE, which may have become quite hot, falling into a PATIENT body cavity.

Subclause 201.15.101 – NEUTRAL ELECTRODES

For low-powered HF SURGICAL EQUIPMENT, for example for dental use, experience has shown that an arrangement where the neutral end of the output circuit is referenced to earth is satisfactory. The return of the HF current from the PATIENT is accomplished capacitively, for example to the earthed metal frame of the dental chair. Consequently this HF SURGICAL EQUIPMENT is exempted from the requirement for a NEUTRAL ELECTRODE.

Subclause 201.15.101.2 – NE cord attachment

The electrical connection of the NE CORD to the part of an NE, except for a MONITORING NE, which is in contact with the PATIENT should be formed such that the NE CONTINUITY MONITOR is capable of detecting any interruption of that connection. MONITORING NES are exempted, since such an interruption is expected to appear as a loss of contact area with the PATIENT.

The test method is suitable for detecting connections which may fuse open during NORMAL USE, however that use is not expected to exceed 1 A.

Subclause 201.15.101.3 – NE cord connector, no conductive parts on PATIENT

In the case of detachment of the NE cord from the NE, it should not be possible for monitoring current from an NE CONTINUITY MONITOR or a CONTACT QUALITY MONITOR to pass through the PATIENT, thus producing a false indication of proper NE attachment.

Subclause 201.15.101.4 – NE cord insulation

Although the voltage difference between the NE application site on the PATIENT and the NE cord conductors may be small, a significant voltage gradient may develop along the PATIENT's body proximal to the surgical site, especially during application of high HF surgical current. Thus, there is a risk of a burn should the NE cord come in contact with a more proximal part of the PATIENT. Application of the HF LEAKAGE CURRENT requirements of 201.8.8.3.102 mitigates this risk. Since lower voltages are expected to be present, the higher LEAKAGE CURRENT limit is deemed appropriate.

Dielectric breakdown of NE cord insulation presents a similar risk to both the PATIENT and the OPERATOR, and thus the HF and mains frequency dielectric strength requirements are deemed necessary. The test voltage magnitudes are unchanged from the prior editions of this particular standard.

Customer: Bradley Hyde- No. of User(s): 1 - Company: Latham and Watkins

Order No.: WS-2023-008299 - IMPORTANT: This file is copyright of IEC, Geneva, Switzerland. All rights reserved.

This file is subject to a licence agreement. Enquiries to Email: sales@iec.ch - Tel.: +41 22 919 02 11

The cables of NEUTRAL ELECTRODES are allowed twice the leakage of the cables of ACTIVE ACCESSORIES, because the voltage levels developed between the conductors of such cables and the PATIENT's skin are generally much lower.

The alternative measured HF leakage capacitance test method may prove simpler to implement than the precedent HF LEAKAGE CURRENT method. See subclause 201.8.3.102 for rationale.

Subclause 201.15.101.5 – NE thermal performance

The literature listed as items [1] through [4] in the bibliography are recommended as a guide in evaluating suitable surrogate surfaces.

This requirement was adopted from ANSI/AAMI HF18:2001, subclause 4.2.3.1. The rationale for that requirement is also adopted, with minor lexical and subclause reference changes for this particular standard, as follows:

The purpose of the NEUTRAL ELECTRODE (NE) in MONOPOLAR electrosurgical procedures is to reliably conduct the required HF surgical current with minimal rise in skin temperature.

Measurements with heated metallic blocks (Moritz & Henriques, 1947 [10]) and with small circular electrodes carrying HF surgical current (Pearce et al., 1983 [12]) show that the maximum safe skin temperature for short-term and long-term exposure is 45 °C. Furthermore, reference to CENELEC Guide 29 [15], Table A1, and interpolating between 48 °C and 43 °C for 8 h or more, gives a maximum allowed surface temperature of 45 °C for 100 min. Normal resting skin temperature varies between about 29 °C and 33 °C, depending on room temperature and humidity. Therefore, NEs that create temperature increases approaching 12 °C cannot be considered safe. Six degrees centigrade represents a conservative safety factor of two and a maximum allowable temperature rise for an acceptable NE. No acceptable NE should exceed a 6 °C temperature rise when subjected to the required current and duration test.

It is recognized that the use of human subjects for qualifying NEs to the requirements of this particular standard may be troublesome or prohibited in many laboratories. However, the specified conformance test is based upon a large volume of empirical data from human tests, using 10 µm infrared imaging instruments, collected and validated by numerous manufacturers and test houses since 1980. Although the use of media and apparatus which yield equivalent results is permitted, documentation of that equivalency must be in place. Therefore, the worst case electrical and thermal properties of NE application sites on a variety of human subjects are the reference standards against which the accuracy of surrogate media and other alternative temperature rise test apparatus are qualified.

Because NE site burns may be confined to very small areas, the qualification measurement must have an adequate spatial sampling frequency to ensure that unacceptable NEs will always be detected. The requirement for one sample per square centimetre is a minimum. Current technology provides for many more samples per square centimetre. However, because noise in the thermal detector can cause individual pixels to appear superheated, a statistical averaging technique should be used to determine the temperature rise within any single square centimetre area. The initial temperature of NEs applied to human skin must be the same in all tests so that all results will be comparable.

At the end of the 60 s application of HF current, the NE is removed from the test surface prior to measuring the final temperature.

HF surgical currents are normally delivered in repetitive short bursts of varying amplitude and duration. Maximum currents and duration of activation depend on the individual technique used and on the type of surgical procedure. The conformance test current is intended to simulate the worst case single activation, with a considerable safety factor. Two sources of information were used to estimate the likely current and duration maxima:

Customer: Bradley Hyde- No. of User(s): 1 - Company: Latham and Watkins

Order No.: WS-2023-008299 - IMPORTANT: This file is copyright of IEC, Geneva, Switzerland. All rights reserved.

This file is subject to a licence agreement. Enquiries to Email: sales@iec.ch - Tel.: +41 22 919 02 11

- A 1973 article in Health Devices presented data in terms of the average currents, voltages, impedances, and minute duty cycles over all procedures studied (ECRI, 1973).
- The unpublished data of Milligan and associates were presented in terms of the maximum, minimum, and average currents and durations for each procedure studied.

These data can be used to estimate population variations. In both studies, it was found that the highest currents and longest durations were found in transurethral (TUR) procedures. For TUR procedures, the ECRI study showed an average CUTTING current of 680 mA and 480 mA for COAGULATION, with duty cycles of 15 % average and 45 % maximum. Milligan studied a smaller sample of 25 TUR procedures performed by 13 surgeons using five electrosurgical units at eight hospitals.

The reported data for all TUR procedures are summarized in the following table. Means and standard deviations σ are calculated over the 25 cases. These data provide useful estimates of the means and variance in measured currents and durations.

Table AA.1 – Summary of measured current and durations for 25 TUR procedures

	Mean	Standard deviation
Length of surgery (hours)	0,86	0,49
Number of activations (per hour)	225	105
CUTTING current		
Maximum current (mA)	407	297
Average current (mA)	297	200
Maximum duration (s)	3,8	2,3
Average duration (s)	2,1	0,7
COAGULATION current		
Maximum current (mA)	339	130
Average current (mA)	258	88
Maximum duration (s)	5,7	7,6
Average duration (s)	2,0	0,7

The total energy dissipated at the NE application site is given by:

$$E = (I_{rms})^2 \times R \times t \text{ where}$$

E is energy dissipated [joules];

I is the NE current [amperes];

t is the duration of current flow [seconds];

R is the real part of the impedance at the NE site [Ω].

The impedance, R , is not generally definable, since its value depends on the NE design and the anatomical structure of the tissue to which it is applied. A "heating factor" Θ may be defined to describe the "stress" placed on an NE as:

$$\Theta = I^2 \times t \text{ (A}^2\text{s)}.$$

Customer: Bradley Hyde- No. of User(s): 1 - Company: Latham and Watkins

Order No.: WS-2023-008299 - IMPORTANT: This file is copyright of IEC, Geneva, Switzerland. All rights reserved.

This file is subject to a licence agreement. Enquiries to Email: sales@iec.ch - Tel.: +41 22 919 02 11

This heating factor has the significance of energy dissipated per Ω of impedance. NES should be able to handle Θ values representative of surgical procedures. A current of 700 mA applied for 60 s yields $\Theta = 30 \text{ A}^2\text{s}$. This value is far in excess of the maximum likely current and duration for a TUR procedure. The maximum likely Θ value can be found by multiplying the square of the largest likely current, i.e. 0,68A from ECRI (1973) [7] data (average) plus one standard deviation, i.e., 0,2A from the Milligan data by the maximum likely duration, i.e., 5,0 s (average) plus one standard deviation, i.e., 7,6 s from the Milligan data, to get

$$\Theta = 9,8 \text{ A}^2\text{s}$$

Thus, 30 A^2s is a conservative test criterion.

A similarly conservative test criterion can be derived for NES marked for "INFANT" use. Since TUR procedures are not performed on infants, a reasonable approach is to use the current and duration data available for general surgical procedures. These data, reported by Pearce (1981), are given in the following table:

Table AA.2 – Summary of measured currents and durations for general surgical procedures

	Mean	Standard deviation
Length of surgery (hours)	1,56	0,84
Number of activations (per hour)	63	84
CUTTING current		
Maximum current (mA)	340	101
Average current (mA)	281	147
Maximum duration (s)	7,6	11
Average duration (s)	2,2	1,8
COAGULATION current		
Maximum current (mA)	267	157
Average current (mA)	198	114
Maximum duration (s)	11	7,5
Average duration (s)	6,5	5,2

Using the data for general surgery and multiplying the square of the maximum likely current plus one standard deviation by the maximum likely duration plus one standard deviation yields

$$\Theta = 3,6 \text{ A}^2\text{s}$$

Thus,

$$\Theta = 15 \text{ A}^2\text{s}$$

is a conservative test criterion and is readily obtained using a current of 500 mA applied for 60 s.

The safety margins inherent in these Θ values are intended to maintain a reasonable margin of safety even in the event of unintended partial loss of contact area between the NE and the PATIENT's skin. Where NES other than MONITORING NES are used, advice to the OPERATOR according to 201.7.9.2.2.101 d) is relied upon to prevent a hazardous loss of contact area. However, where CONTACT QUALITY MONITORS and MONITORING NES are in use, the OPERATOR expects to be relieved of the burden of NE contact surveillance, relying fully upon the CONTACT QUALITY MONITOR to alert the OPERATOR to area loss before it becomes hazardous.

Customer: Bradley Hyde- No. of User(s): 1 - Company: Latham and Watkins

Order No.: WS-2023-008299 - IMPORTANT: This file is copyright of IEC, Geneva, Switzerland. All rights reserved.

This file is subject to a licence agreement. Enquiries to Email: sales@iec.ch - Tel.: +41 22 919 02 11

Therefore, MONITORING NES are tested with the same area loss which will cause the CONTACT QUALITY MONITOR to sound an alarm.

References are found in the bibliography items [7] through [12].

The test currents by weight range in Table 201.3 were derived as follows. NEUTRAL ELECTRODES for adults, when tested with the 700 mA current based on the HF18 standard, result in a heating factor of 30 A²s.

NEUTRAL ELECTRODES for children (PATIENT weight 5 kg to 15 kg) have active contact areas that are approximately one half of adult sized NEUTRAL ELECTRODES. When tested with the 500 mA current based on the HF18 standard, results in a heating factor of 15 A²s, which is one half of the maximum allowed adult value.

NEUTRAL ELECTRODES for newborns (PATIENT weight less than 5 kg) have active contact areas that are approximately one half as those for children and thus using a heating factor that is one half of that used for children was chosen. This results in a heating factor of 7,5 A²s which implies the test current of 350 mA. Although there is no statistical data to prove the selection of this test current, surgical power settings for these small PATIENTS are always very low so it is felt the 350 mA test current for 60 s results in a reasonable safety margin.

Subclause 201.15.101.6 – NE contact impedance

This requirement was adopted from ANSI/AAMI HF18:2001, subclause 4.2.3.2. The 200 kHz phase angle criterion for distinguishing conductive and capacitive NEs was developed, lacking a clear published definition, *a priori* or otherwise.

The rationale from ANSI/AAMI HF18:2001, subclause A.4.2.3.2 is also adopted, with minor lexical and subclause reference changes for this particular standard, as follows:

The contact impedance must be low enough that the NEUTRAL ELECTRODE represents the preferred current pathway. In the case of HF SURGICAL EQUIPMENT having an EARTH REFERENCED PATIENT CIRCUIT, this will minimize the possibility of alternate return current paths other than via the NE. A value of 75 Ω is judged an acceptable maximum contact impedance for conductive NES when measured according to ANSI/AAMI HF18:2001 using human subjects. However, that standard imposes a 50 Ω limit when a metal plate is used in lieu of a human subject; this reduction compensates for the impedance contribution of the deeper subcutaneous tissue which becomes part of the measured NE contact impedance.

Since the impedance of capacitive NES varies as the inverse of the frequency, it is appropriate to describe their impedance characteristics in terms of capacitance. A value of 4 nF was specified as the minimum acceptable capacitance because it is consistent with the characteristics of the majority of capacitive NES which have been commercially available for many years and found to be clinically acceptable.

The test current of 200 mA represents the low limit of average currents from the two studies cited above. Tissue-NE impedance generally increases as the current decreases, making the lower limit preferable. The frequency range of 200 kHz to 5 MHz is believed to encompass the range over which MONOPOLAR HF SURGICAL EQUIPMENT develops significant energy levels.

The dimensions of the metallic test plate should be at least as large as the NEUTRAL ELECTRODE.

Capacitive NES are permitted a higher impedance because they do not dissipate heat.

Subclause 201.15.101.7 – NE adhesion

This requirement was adopted from ANSI/AAMI HF18:2001, subclause 4.2.3.3.

After application, NES, except MONITORING NES, should remain in place when subjected to stresses that may occur during customary use as a result of the site chosen for placement, inadvertent pulling, or accidental contact with preparatory solutions or physiologic fluids. MONITORING NES are exempt from this requirement because contact area loss due to adhesive failure is expected to cause a CONTACT QUALITY MONITOR alarm, thus preventing a hazard to the PATIENT.

Subclause 201.15.101.8 – NE shelf life

The adhesives and conductive gels used on single use NES may deteriorate over time, even when stored according to the instructions for use. Therefore, it is necessary to determine that these devices conform after storage until the marked expiration date.

Clause 202

HF surgery is a very long established modality with known interference inherent during activation. Since the clinical benefits of HF SURGICAL EQUIPMENT outweigh the risks of interference and since HF SURGICAL EQUIPMENT is normally operated for short periods only, this type of equipment is exempt from the requirements of IEC 60601-1-2 subclause 6.1.1.1 when it is activated.

HF SURGICAL EQUIPMENT performs its CUTTING and COAGULATION functions through the use of radio frequency energy, and HF EMISSION frequently much above the CISPR 11 limits is present. The power levels and harmonic content of the output of the HF SURGICAL EQUIPMENT are necessary to enable the HF SURGICAL EQUIPMENT to carry out its clinical function effectively.

The EMISSIONS strongly depend on the arrangement and length of the active and neutral cords, on the operating mode (sparking or not) and on many other application conditions. Furthermore, many diagnostic, monitoring, anaesthetic and infusion EQUIPMENT have APPLIED PARTS or PATIENT circuits which are directly connected to the PATIENT. For such equipment, particular test arrangements simulating direct connection to the PATIENT circuits of a HF SURGICAL EQUIPMENT are necessary for testing EMC IMMUNITY (see 202.6.2.101 and Figures 201.109, 201.110 and 201.111 of IEC 60601-2-34⁴⁾). This was considered the best way to assure electromagnetic compatibility between HF SURGICAL EQUIPMENT and some other medical devices used in its vicinity.

For a standardized source of electromagnetic interference to be used for such tests, the following conditions have been determined in IEC 60601-2-34:

"The HF SURGICAL EQUIPMENT shall comply with IEC 60601-2-2, shall have a minimum cut mode of 300 W, a minimum COAGULATION mode of 100 W and a working frequency of 450 kHz $\pm$ 100 kHz."

However, during stand-by operation, the HF SURGICAL EQUIPMENT may be energized for long periods and compliance with the EMC requirements is considered necessary.

During the immunity tests of IEC 61000-4-3 and IEC 61000-4-6, the manufacturer will need to specify how compliance to the standard is checked. This includes precautions needed to

⁴⁾ Third edition, in preparation: IEC 60601-2-34, *Medical electrical equipment – Part 2-34: Particular requirements for the basic safety and essential performance of invasive blood pressure monitoring equipment*

ensure the duty cycle of the generator is not exceeded as well as how perturbations in the output power are detected.

Additional information on the electromagnetic EMISSIONS created by HF SURGICAL EQUIPMENT may be found in Annex BB.

Withdrawn

Annex BB (informative)

ELECTROMAGNETIC DISTURBANCES created by HF SURGICAL EQUIPMENT

BB.1 Overview

Medical devices used in surgery are exposed to many sources of ELECTROMAGNETIC DISTURBANCE (EMD). The most prevalent source is from the HF SURGICAL EQUIPMENT that is used to cut and coagulate tissue. Although there are standards for many types of EMD, there is little information available regarding the ELECTROMAGNETIC DISTURBANCE created by HF SURGICAL EQUIPMENT.

The purpose of this annex is to provide medical device manufacturers with information about the specific types and levels of EMISSIONS generated by HF SURGICAL EQUIPMENT. It also includes tests which manufacturers may wish to use in determining if their designs are resistant to these type of ELECTROMAGNETIC DISTURBANCE.

BB.2 Terms and definitions

For the purposes of this annex, the definitions of terms appearing in small capitals are those from this particular standard and the standards listed in clause 1.3 of the general standard, plus the following terms and definitions.

NOTE Definitions for ELECTROMAGNETIC DISTURBANCE and EMISSIONS may be found in IEC 60601-1-2.

BB.2.1

E-FIELD

electric field present in the far field as induced by the magnetic field from HF SURGICAL EQUIPMENT

BB.2.2

H-FIELD

magnetic field induced by the flow of current from the HF SURGICAL EQUIPMENT

BB.3 Technical information

BB.3.1 General information about HF SURGICAL EQUIPMENT

During surgery, HF energy may be used for CUTTING tissue or to provide hemostasis (COAGULATION). This energy is generated by the HF SURGICAL EQUIPMENT and delivered to the surgical site using various sterile ACCESSORIES. The frequency of the HF energy is typically between 200 kHz and 1 MHz. These frequencies are high enough that human tissue cannot respond to them, and thus no nerve or muscle stimulation occurs. All of the surgical effect is due to the current density of the HF energy.

The HF energy may be delivered to the surgical site in one of two ways. The first method is called MONOPOLAR or unipolar. This means that the surgical effect occurs at a single pole which is under the surgeon's control. The energy is generated in the HF SURGICAL EQUIPMENT, is carried through a cord to an ACCESSORY held by the surgeon, through the PATIENT, is collected by a large surface area PATIENT return electrode (NEUTRAL ELECTRODE) and is carried back to the HF SURGICAL EQUIPMENT. It is the current density at the tip of the ACCESSORY ACTIVE ELECTRODE(s) that causes the localized surgical effect. After entering the PATIENT's body, the current disperses, limiting the area of the surgical effect. The large surface area of the NEUTRAL ELECTRODE is designed to keep the current density low to prevent heating or other tissue effects. The PATIENT return electrode is the second pole in the circuit. The most

Customer: Bradley Hyde- No. of User(s): 1 - Company: Latham and Watkins

Order No.: WS-2023-008299 - IMPORTANT: This file is copyright of IEC, Geneva, Switzerland. All rights reserved.

This file is subject to a licence agreement. Enquiries to Email: sales@iec.ch - Tel.: +41 22 919 02 11

SERENDIA-00084220

Ex. 2077-072

CX-2848.72

common MONOPOLAR ACCESSORY is the HF surgical pencil, so named because it resembles a thick pencil held by the surgeon.

The second method of energy delivery is called BIPOLAR. The surgical ACCESSORY used by the surgeon has two electrodes, each with a small surface area. The HF energy passes from the HF surgical unit to one electrode, through the tissue, to the other electrode and back to the HF surgical unit. The area of the electrodes and the tissue between them is small and so the current density is high. Thus, the surgical effect occurs only in the tissue grasped between the electrodes. A NEUTRAL ELECTRODE is not required. The most common BIPOLAR ACCESSORY is HF surgical forceps.

Most HF SURGICAL EQUIPMENT allows the user to control the output power as a means of controlling the depth and speed of the surgical effect. The output voltage and current may vary depending on the HF SURGICAL MODE, the power setting and the load presented to the HF SURGICAL EQUIPMENT.

The surgical effect of CUTTING is generally achieved using a sine wave with a voltage between 200 V and 1 200 V. The current density at the tip of the electrode causes heating of the contents of cells immediately adjacent to the electrode. The cell contents turn to steam and the cell wall ruptures. The electrode moves through this steam layer and very small arcs pass from the electrode tip to the tissue. A pure sine wave cuts with little or no hemostasis. If the sine wave is interrupted, various levels of hemostasis may be achieved in addition to the CUTTING action. The lower the duty cycle the greater the hemostasis. However, lowering the duty cycle also requires that the peak voltage be increased to achieve the same output power. Power levels used in the cut mode range between 10 W and 300 W.

The surgical effect of COAGULATION may be achieved using several different methods. A pure sine wave which is below 200 V will not cut tissue but will desiccate and coagulate tissue. This waveform does not produce arcs. It is used for contact COAGULATION in both the MONOPOLAR and BIPOLAR modes. When the surgeon needs to coagulate bleeding tissue without touching it, a high voltage interrupted sinusoidal waveform is generally used. This waveform may use a voltage between 1 200 V and 4 600 V. Power levels used for the MONOPOLAR COAGULATION mode range from 10 W to 120 W. Power levels for the BIPOLAR COAGULATION mode range from 1 W to 100 W.

The worst case EMISSIONS created by HF SURGICAL EQUIPMENT occur during activation of the COAGULATION mode at the maximum power setting while arcing to tissue or metal.

BB.3.2 Types of EMISSIONS created by HF SURGICAL EQUIPMENT

BB.3.2.1 Radiated

During surgery, the therapeutic current flows from the HF surgical unit through the ACCESSORY cable, through the PATIENT, through an ACCESSORY cable again, and back to the unit. This circuit may take on different forms, sizes and arrangements. The current flowing creates both a radiated E-FIELD and an H-FIELD. These fields may couple to the ACCESSORY, or POWER SUPPLY CORD used by other equipment. The worst case scenario for E-FIELD coupling is to have the HF SURGICAL ACCESSORY cables in close proximity and parallel with other ACCESSORY cables. E-FIELD coupling is also made worse during clinical situations where arcs occur. The worst case scenario for H-FIELD coupling is to have the HF surgical circuit spread out in a large circle and other ACCESSORY cables attached to the PATIENT who is within that circle. E-FIELD coupling typically generates worst case EMISSIONS that are higher in frequency (tens to hundreds of megahertz) than H-FIELD coupling (tens to hundreds of kilohertz).

BB.3.2.2 Conducted through the mains POWER SUPPLY CORD

Electromagnetic noise conducted through the mains POWER SUPPLY CORD increases during activation of the HF SURGICAL EQUIPMENT through a combination of internal coupling to the HF output and high voltage power supplies that are only active during HF output generation.

Customer: Bradley Hyde- No. of User(s): 1 - Company: Latham and Watkins

Order No.: WS-2023-008299 - IMPORTANT: This file is copyright of IEC, Geneva, Switzerland. All rights reserved.

This file is subject to a licence agreement. Enquiries to Email: sales@iec.ch - Tel.: +41 22 919 02 11

BB.3.2.3 Conducted through the PATIENT

The therapeutic current that is applied to the PATIENT to achieve CUTTING and COAGULATION impresses a voltage on the PATIENT that may be coupled into other equipment. This coupling may be direct or capacitive. Direct coupling occurs into the inputs of devices that are measuring PATIENT voltages (e.g. ECG, EEG, EMG, evoked potential monitors). Capacitive coupling occurs when equipment cables or sensors are in close contact with the PATIENT (e.g. pulse oximeter probes, invasive blood pressure transducers, temperature probes, camera systems). A combination of these methods is possible. The value of the voltage impressed on the PATIENT is highly dependant on the HF SURGICAL MODE used. BIPOLAR modes utilize peak-to-peak voltages ranging from tens to a few hundred of volts and generate little or no sparking. CUTTING modes utilize peak to peak voltages from several hundred to a few thousand volts and generate very small sparks. COAGULATION modes utilize peak to peak voltages from a few thousand up to fourteen thousand volts with large sparks frequently being desired. Generally only a fraction of the HF voltage is coupled into other equipment but for devices that measure in the millivolt or microvolt range, that can be a problem.

BB.3.3 Measurement techniques


For the purpose of this annex, the measurements were taken using techniques intended to create the worst case values that may be experienced by MEDICAL ELECTRICAL EQUIPMENT during surgery.

The measurements reported below were taken multiple times using all of the output modes available and using the maximum output powers the units were capable of. Four different clinical situations were simulated. These situations were: open circuit activation, activation at the RATED LOAD of the HF SURGICAL EQUIPMENT (the load which produces the maximum output power), sparking to metal, and sparking to a saline-soaked sponge to simulate sparking to tissue.

All of these measurements were repeated multiple times using HF SURGICAL EQUIPMENT from a variety of manufacturers. The resulting data were used to create the worst case values of BB.3.4.4.

BB.3.3.1 E-FIELD measurements

A non-conductive table 1 m above a ground plane was used to support the ACCESSORY cables from the HF SURGICAL EQUIPMENT under test. The measurement techniques found in CISPR 11 were used. The set up is illustrated in Figure BB.1. The measurements were recorded as peak or quasi-peak values that occur between 30 MHz and 1 GHz.

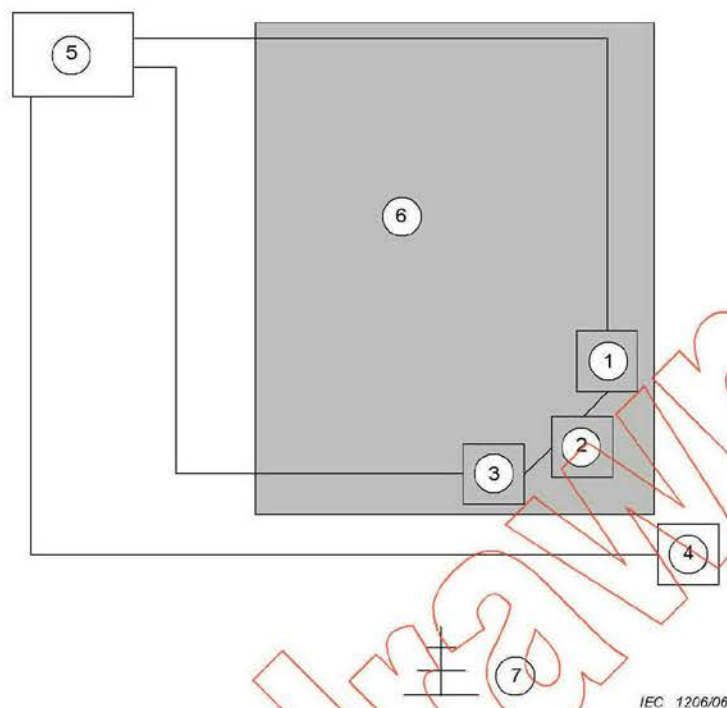


- Key**
- ① ACTIVE ACCESSORY
 - ② Load
 - ③ NE or saline soaked sponge
 - ④ Footswitch
 - ⑤ HF SURGICAL EQUIPMENT
 - ⑥ Non-conductive table
 - ⑦ Antenna – 10 m distance, vertical polarity

Figure BB.1 – E-FIELD EMISSIONS test setup

BB.3.3.2 H-FIELD measurements

A non-conductive table 1 m above a ground plane was used to support the ACCESSORY cables from the HF SURGICAL EQUIPMENT unit under test. The set up is illustrated in Figure BB.2.

The measurements were recorded as peak or quasi-peak values that occur between 10 kHz and 30 MHz.

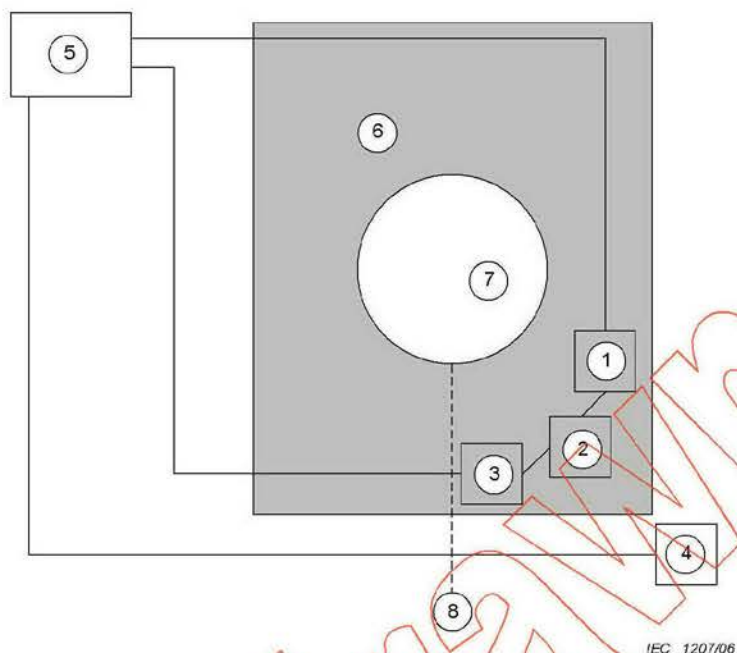

- Key**
- ① ACTIVE ACCESSORY
 - ② Load
 - ③ NE or saline soaked sponge
 - ④ Footswitch
 - ⑤ HF SURGICAL EQUIPMENT
 - ⑥ Non-conductive table
 - ⑦ Antenna
 - ⑧ Cable to measurement equipment

Figure BB.2 – H-FIELD EMISSIONS test setup

BB.3.3.3 Mains conducted measurements

A non-conductive table 1 m above a ground plane was used to support the ACCESSORY cables from the HF SURGICAL EQUIPMENT under test. The set up is illustrated in Figure BB.3.

The measurements were recorded as peak or quasi-peak values that occur between 150 kHz and 30 MHz.

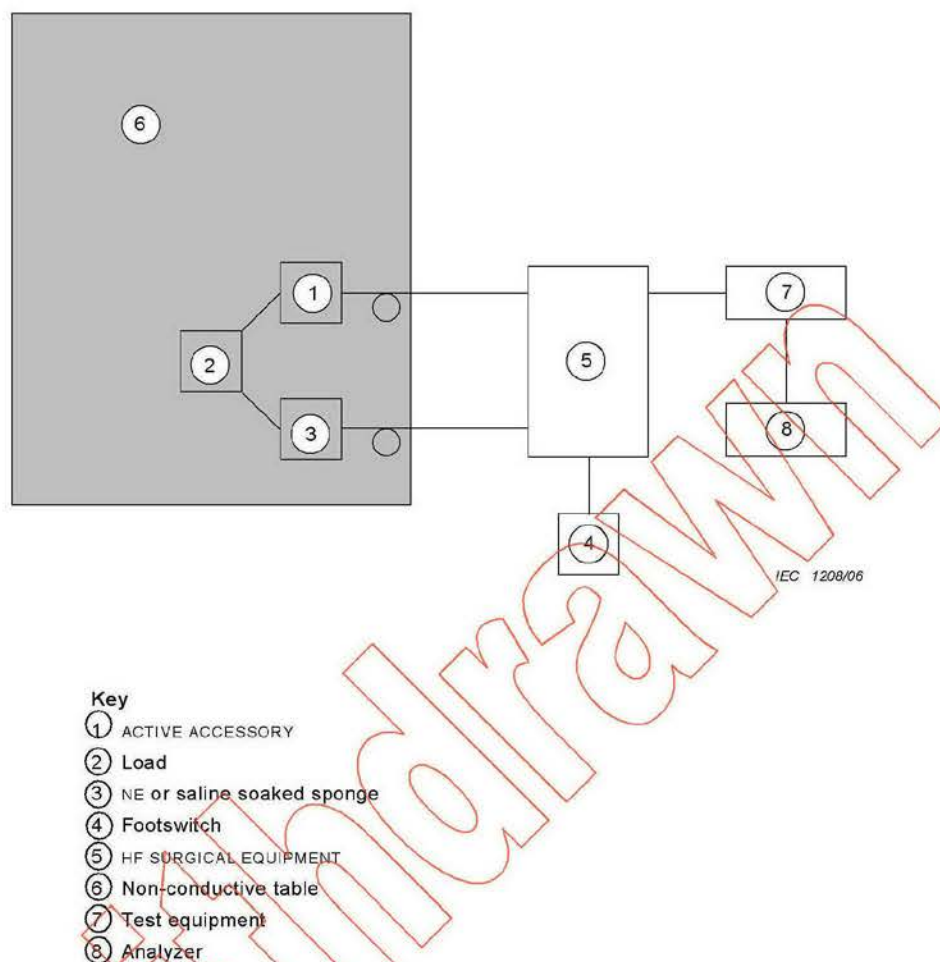


Figure BB.3 – Conducted EMISSIONS test setup

BB.3.4 Data summary**BB.3.4.1 E-FIELD EMISSIONS**

The greatest values were typically below 50 MHz, with lower energy at higher frequencies. Arcing increases the energy at all frequencies, with arcing to metal being the worst case clinical situation.

BB.3.4.2 H-FIELD EMISSIONS

The greatest values were typically at the fundamental frequency of the HF SURGICAL EQUIPMENT, with additional peaks at multiples of the fundamental frequency. Arcing increases the energy at all frequencies, with arcing to metal being the worst case clinical situation.

BB.3.4.3 Mains conducted EMISSIONS

The greatest values were typically at the fundamental frequency of the HF SURGICAL EQUIPMENT with additional peaks at multiples of the fundamental frequency. Arcing increases the energy at all frequencies, with arcing to metal being the worst case clinical situation.

Customer: Bradley Hyde- No. of User(s): 1 - Company: Latham and Watkins

Order No.: WS-2023-008299 - IMPORTANT: This file is copyright of IEC, Geneva, Switzerland. All rights reserved.

This file is subject to a licence agreement. Enquiries to Email: sales@iec.ch - Tel.: +41 22 919 02 11

BB.3.4.4 Maximum EMISSION levels of HF SURGICAL EQUIPMENT

The greatest level of EMISSIONS was generated by spark gap units. This type of HF SURGICAL EQUIPMENT is no longer sold but still found in many hospitals. This type of unit creates the worst case EMD environment due to a very high output voltage and the use of a spark gap to create COAGULATION waveforms. The use of a spark gap tends to generate much higher levels of EMISSIONS at higher frequencies. The worst case EMISSION values are shown in Table BB.1 and Table BB.2. The E-FIELD measurements were from a distance of 10 m.

Table BB.1 – Worst case emissions of spark gap type HF SURGICAL EQUIPMENT

EMISSION type	No arcing	Arcing to saline	Arcing to metal
E-FIELD	92 dB μ V/m (40 mV/m)	80 dB μ V/m (10 mV/m)	95 dB μ V/m (56 mV/m)
H-FIELD	96,47 dB μ A/m (67 mA/m)	99,47 dB μ A/m (94 mA/m)	96,47 dB μ A/m (67 mA/m)
Mains conducted	117 dB μ V (708 mV)	Not measured	Not measured

Table BB.2 – Worst case emissions of non-spark gap (modern) HF SURGICAL EQUIPMENT

EMISSION type	No arcing	Arcing to saline	Arcing to metal
E-FIELD	78 dB μ V/m (8mV/m)	77 dB μ V/m (7 mV/m)	83 dB μ V/m (14 mV/m)
H-FIELD	61,47 dB μ A/m (1,1 mA/m)	63,47 dB μ A/m (1,5 mA/m)	62,47 dB μ A/m (1,3 mA/m)
Mains conducted	97 dB μ V (71 mV)	Not measured	100 dB μ V (100 mV)


BB.4 Suggested tests

BB.4.1

The following information describes some ad hoc tests that have been used by equipment manufacturers to determine if their products can withstand the EMISSIONS produced by HF SURGICAL EQUIPMENT. These tests are meant to serve as guides only and may be modified based on how the equipment is situated with respect to the HF SURGICAL EQUIPMENT. The tests below were designed to simulate the two types of equipment being situated in close proximity (both the ENCLOSURES and the cables). Just as in IEC 60601-1-2, the equipment manufacturer should define what the acceptable response to this test should be prior to conducting it.

BB.4.2

Set up the equipment that is to be tested. Wrap the cord of a MONOPOLAR HF SURGICAL ACCESSORY around the equipment so that at least two full loops of the cord are present as shown in Figure BB.4.

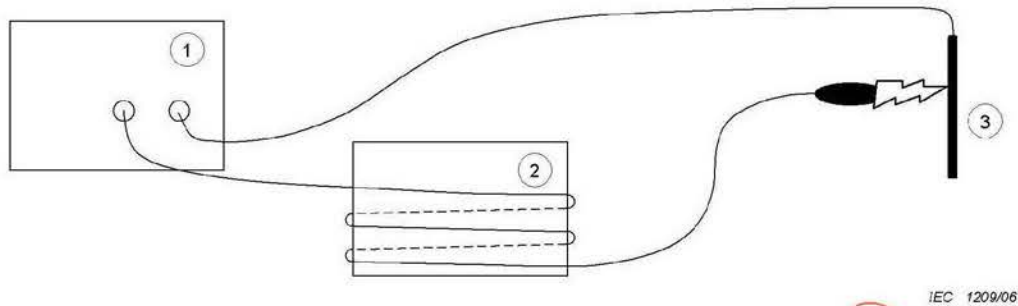**Key**

- ① HF SURGICAL EQUIPMENT
- ② Unit under test
- ③ Metal plate

Figure BB.4 – Unit ad hoc test

Attach one end of a cord to the NEUTRAL ELECTRODE connector of the HF SURGICAL EQUIPMENT and the other end to a metal plate. Using the MONOPOLAR HF SURGICAL ACCESSORY, activate the HF SURGICAL EQUIPMENT in each possible output mode and arc the ACCESSORY to the metal plate. For each mode, adjust the HF SURGICAL EQUIPMENT to the setting which will create the highest peak output voltage.

This test generates high E-FIELDS and high H-FIELDS with the greatest possible spread of frequencies.

BB.4.3

Repeat the test of BB.4.2 with the MONOPOLAR HF SURGICAL ACCESSORY short circuited to (touching) the metal plate. The HF SURGICAL EQUIPMENT should be adjusted to obtain the maximum output power for each output mode.

This test generates the highest output currents and thus the greatest H-FIELDS. It also creates high E-FIELDS at the fundamental output frequency.

BB.4.4

Repeat the tests of BB.4.2 and BB.4.3 with the cord of the MONOPOLAR HF SURGICAL ACCESSORY wrapped around the mains power cord of the unit under test as shown in Figure BB.5.

This test simulates the noise that can be coupled into the equipment through the mains power cord.

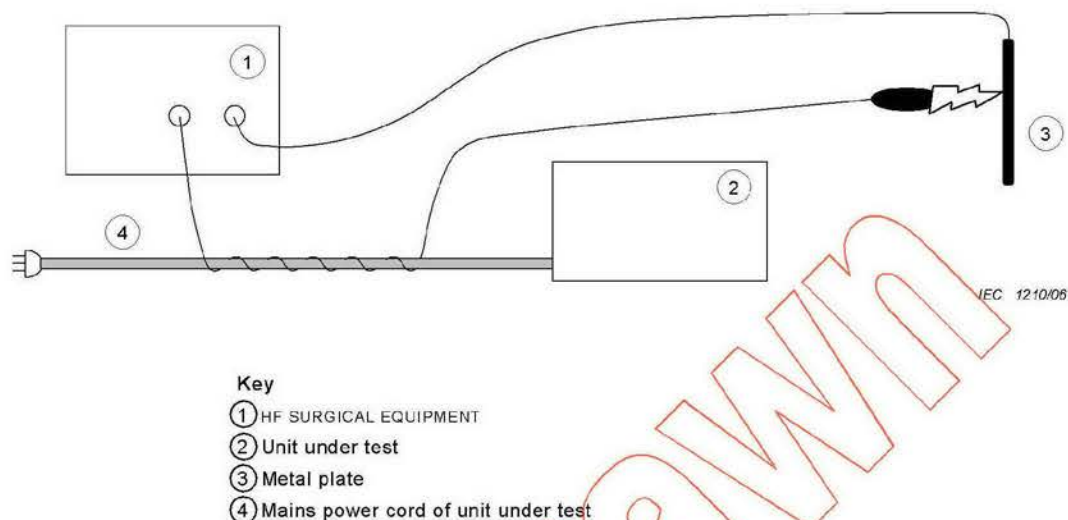


Figure BB.5 – Power cord ad hoc test

BB.4.5

If the equipment has cords that enter the sterile field, coupling can also occur between those cords and the MONOPOLAR HF SURGICAL ACCESSORY cord. To test for this possibility, repeat the tests of BB.4.2 and BB.4.3 with the cord of the MONOPOLAR HF SURGICAL ACCESSORY wrapped around the ACCESSORY cord of the unit under test as shown in Figure BB.6.

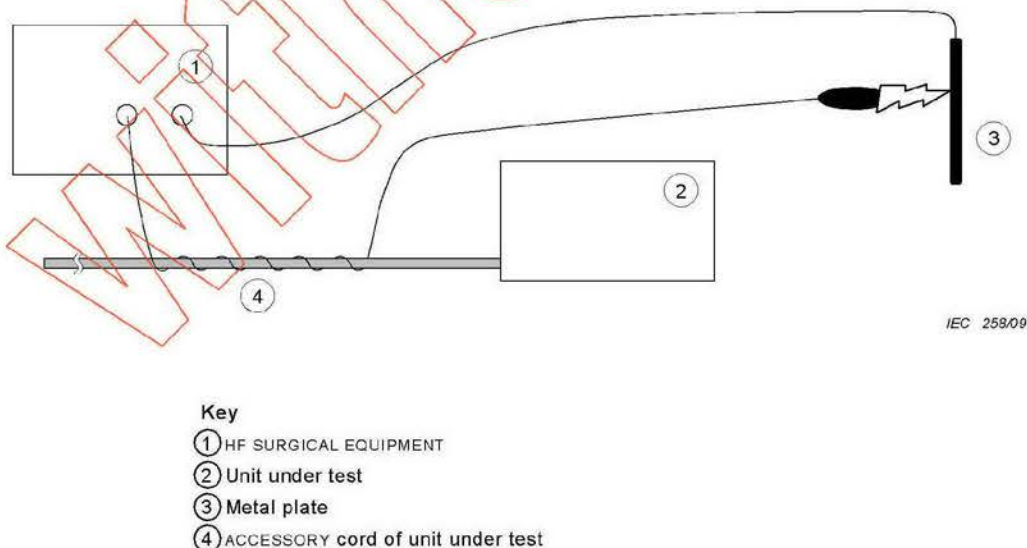


Figure BB.6 – ACCESSORY cord ad hoc test

BB.4.6

Tests to determine the impact of EMISSIONS that are conducted through the PATIENT may vary widely based on how well coupled the equipment is to the PATIENT. The reader is urged to consult the particular standard(s) for their type of equipment for additional information. Many of these particular standards have already included this type of test.

Customer: Bradley Hyde- No. of User(s): 1 - Company: Latham and Watkins

Order No.: WS-2023-008299 - IMPORTANT: This file is copyright of IEC, Geneva, Switzerland. All rights reserved.

This file is subject to a licence agreement. Enquiries to Email: sales@iec.ch - Tel.: +41 22 919 02 11

Bibliography

- [1] NESSLER N., REISCHER W., SALCHNER M. *Measurement Science Review*, 2003, Volume 3, Section 2.
- [2] NESSLER N. *Current Density distribution in Human skin under the Grounding electrode of Electrosurgery*, BEMS 17th Annual Meeting, Boston, MA., 1995
- [3] NESSLER N. HUTER H., WANG L. Sicherheitstester für HF-Chirurgie-Neutralelektroden. *Biomedizinische Technik*, 1993, Volume 38, pp 5-9.
- [4] NESSLER N. REISCHER W., SALCHNER M. *Electronic Skin – Test Device For Electrosurgical Electrodes*. 12th IMEKO TC4 International Symposium, Zagreb 2002
- [5] KELLER A., ROSENFELDER G. "DIN EN 60601-2-2, 4th edition, clause 59.103.5 / 59.104.4 – comparison of alternative test methods – leakage current test method versus capacitance test method", Aesculap AG & Co. KG, 25-Aug-05.
- [6] KÖNIG A., HEINRICH M. Comparative test of HF leakage current on cables with different measuring methods. *BOWA Electronic*, 11.03.05.
- [7] EMERGENCY CARE RESEARCH INSTITUTE. Clinical studies. *Health Devices*, 1973, volume 2, numbers 8-9, pp. 194-195.
- [8] EMERGENCY CARE RESEARCH INSTITUTE. *Draft Environmental Requirements and Test Methods for Non-Implantable Medical Devices* (Contract No. FDA-74-230). Plymouth Meeting, PA: ECRI, July 1978.
- [9] EMERGENCY CARE RESEARCH INSTITUTE. *Development of Environmental Test Methods for Non-Implantable Medical Devices, Final Report* (Contract No. 223-77-5035). Plymouth Meeting, PA: ECRI, April 1979.
- [10] MORITZ, AR, HENRIQUES, FC. Studies in thermal injury: II. The relative importance of time and surface temperature in the causation of cutaneous burns. *American Journal of Pathology*, 1947, volume 23, number 5, pp. 695-720.
- [11] PEARCE, JA, FOSTER, KS, MULLIKIN, JC, GEDDES, LA. *Investigations and Studies on Electrosurgery*, (HHS publication FDA 84-4186). Rockville, MD: U.S. Food and Drug Administration, 1981.
- [12] PEARCE, JA, GEDDES, LA, VAN VLEET, JF, FOSTER, K, ALLEN, J. Skin burns from electrosurgical current. *Medical Instrumentation*, 1983, volume 17, number 3, pp. 225-231.
- [13] IEC 60601-2-18:1996, *Medical electrical equipment – Part 2-18: Particular requirements for the safety of endoscopic equipment*
Amendment 1: 2000
- [14] ANSI/AAMI HF18:2001, *Electrosurgical Devices*
- [15] CENELEC Guide 29:2007, *Temperatures of hot surfaces likely to be touched*

SOMMAIRE

AVANT-PROPOS.....	82
INTRODUCTION.....	85
201.1 Domaine d'application, objet et normes connexes.....	86
201.2 Références normatives.....	88
201.3 Termes et définitions.....	88
201.4 Exigences générales.....	92
201.5 Exigences générales relatives aux essais des APPAREILS EM.....	92
201.6 Classification des APPAREILS EM et des SYSTÈMES EM.....	92
201.7 Identification, marquage et documentation DES APPAREILS EM.....	93
201.8 Protection contre les DANGERS d'origine électrique provenant des APPAREILS EM.....	98
201.9 Protection contre les dangers mécaniques des APPAREILS EM et SYSTÈMES EM.....	114
201.10 Protection contre les DANGERS dus aux rayonnements involontaires ou excessifs.....	114
201.11 Protection contre les températures excessives et les autres DANGERS.....	114
201.12 Précision des commandes et des instruments, et protection contre les caractéristiques de sortie présentant des dangers.....	116
201.13 SITUATIONS DANGEREUSES et conditions de défaut.....	122
201.14 SYSTÈMES ÉLECTROMÉDICAUX PROGRAMMABLES (SEMP).....	122
201.15 Construction de l'APPAREIL EM.....	122
201.16 SYSTÈMES EM.....	127
201.17 Compatibilité électromagnétique des APPAREILS EM ET DES SYSTÈMES EM.....	127
202 * Compatibilité électromagnétique – Exigences et essais.....	127
208 Exigences générales, essais et guide pour les systèmes d'alarme des appareils et des systèmes électromédicaux.....	128
Annexes.....	128
Annexe AA (informative) Guide particulier et justifications.....	129
Annexe BB (informative) PERTURBATIONS ÉLECTROMAGNÉTIQUES GÉNÉRÉES PAR LES APPAREILS D'ÉLECTROCHIRURGIE HF.....	153
Bibliographie.....	162
Figure 201.101 – Symbole utilisé avec un CIRCUIT PATIENT RÉFÉRENCÉ À LA TERRE.....	93
Figure 201.102 – Symbole utilisé avec un CIRCUIT PATIENT ISOLÉ HF.....	93
Figure 201.103 – Circuit adapté aux essais de conformité selon 201.8.4.101.....	99
Figure 201.104 – Mesurage du COURANT DE FUITE HF avec ÉLECTRODE NEUTRE référencée par rapport à la terre et résistance de charge entre électrodes.....	102
Figure 201.105 – Mesurage du COURANT DE FUITE HF avec ÉLECTRODE NEUTRE référencée par rapport à la terre et résistance de charge entre l'ÉLECTRODE ACTIVE et la terre.....	103
Figure 201.106 – Mesurage du COURANT DE FUITE HF avec ÉLECTRODE NEUTRE isolée de la terre à HAUTE FRÉQUENCE.....	104
Figure 201.107 – Mesurage du COURANT DE FUITE HF d'une ÉLECTRODE BIPOLAIRE.....	105
Figure 201.108 – Montage d'essai pour dispositifs d'arrêt de traction et de torsion de câbles d'ACCESSOIRE ACTIF.....	112
Figure 201.109 – Mesurage de la puissance de sortie (assignée) – sortie MONOPOLAIRE.....	117

Customer: Bradley Hyde- No. of User(s): 1 - Company: Latham and Watkins

Order No.: WS-2023-008299 - IMPORTANT: This file is copyright of IEC, Geneva, Switzerland. All rights reserved.

This file is subject to a licence agreement. Enquiries to Email: sales@iec.ch - Tel.: +41 22 919 02 11

Figure 201.110 – Mesurage de la puissance de sortie (assignée) – sortie BIPOLAIRE	118
Figure 201.111 – Méthode d'essai de rétroaction d'une sortie active à une autre dans une activation simultanée	121
Figure AA.1 – Exemple de différentes parties d'un APPAREIL D'ÉLECTROCHIRURGIE HF	130
Figure AA.2 – FACTEUR DE CRÊTE par rapport à la tension de crête.....	135
Figure AA.3 – Exemple de circuit PATIENT avec ÉLECTRODE NEUTRE référencée par rapport à la terre aux fréquences de fonctionnement	139
Figure BB.1 – Montage d'essai pour les ÉMISSIONS DE CHAMP-E.....	156
Figure BB.2 – Montage d'essai pour les ÉMISSIONS DE CHAMP-H.....	157
Figure BB.3 – Montage d'essai pour ÉMISSIONS conduites	158
Figure BB.4 – Essai pratique (ad hoc) de l'unité	160
Figure BB.5 – Essai pratique (ad hoc) du câble d'alimentation.....	161
Figure BB.6 – Essai pratique (ad hoc) du câble de l'ACCESSOIRE	161
Tableau 201.101 – Couleurs des voyants lumineux et leur signification pour les APPAREILS D'ÉLECTROCHIRURGIE HF.....	94
Tableau 201.102 – Puissances de sortie maximales dans des CONDITIONS DE PREMIER DÉFAUT	120
Tableau 201.103 – Courants d'essai par gamme de poids	124
Tableau AA.1 – Synthèse des intensités de courant et durées mesurées pour 25 interventions TUR	148
Tableau AA.2 – Synthèse des intensités de courant et durées mesurées pour des interventions chirurgicales générales	149
Tableau BB.1 – Emissions du cas le plus défavorable d'APPAREIL D'ÉLECTROCHIRURGIE HF de type à éclateur	159
Tableau BB.2 – Emissions du cas le plus défavorable d'APPAREIL D'ÉLECTROCHIRURGIE HF de type sans éclateur (moderne)	159

Customer: Bradley Hyde- No. of User(s): 1 - Company: Latham and Watkins

Order No.: WS-2023-008299 - IMPORTANT: This file is copyright of IEC, Geneva, Switzerland. All rights reserved.

This file is subject to a licence agreement. Enquiries to Email: sales@iec.ch - Tel.: +41 22 919 02 11

SERENDIA-00084231

Ex. 2077-083

CX-2848.83