




US 20220145380A1

(19) **United States**(12) **Patent Application Publication**
DOAN et al.(10) **Pub. No.: US 2022/0145380 A1**(43) **Pub. Date: May 12, 2022**(54) **COST-EFFECTIVE DETECTION OF LOW
FREQUENCY GENETIC VARIATION****Publication Classification**(71) Applicant: **CHILDREN'S MEDICAL CENTER
CORPORATION**, Boston, MA (US)(51) **Int. Cl.**
C12Q 1/6851 (2006.01)
C12Q 1/6886 (2006.01)(72) Inventors: **Ryan N. DOAN**, Boston, MA (US);
Christopher A. WALSH, Boston, MA
(US)(52) **U.S. Cl.**
CPC **C12Q 1/6851** (2013.01); **C12Q 2600/156**
(2013.01); **C12Q 2600/118** (2013.01); **C12Q**
1/6886 (2013.01)(73) Assignee: **CHILDREN'S MEDICAL CENTER
CORPORATION**, Boston, MA (US)(21) Appl. No.: **17/427,394**(57) **ABSTRACT**(22) PCT Filed: **Nov. 26, 2019**(86) PCT No.: **PCT/US2019/063447**

§ 371 (c)(1),

(2) Date: **Jul. 30, 2021****Related U.S. Application Data**(60) Provisional application No. 62/799,671, filed on Jan.
31, 2019.

Methods are described for the detection of low frequency genetic variants, such as somatic mosaic variants. The methods comprise parallel amplification reactions of a target nucleic acid sequence to generate overlapping amplicons, pooled sequencing of the amplicons, and demultiplexed detection of low frequency variants.

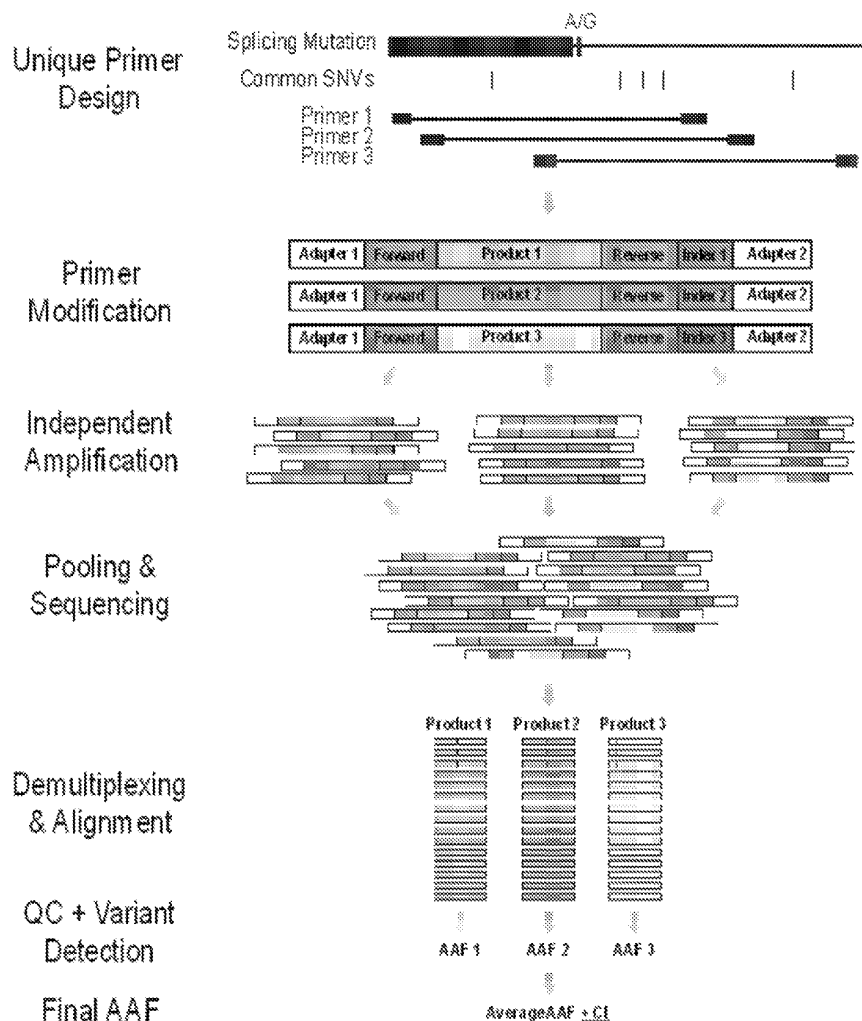Specification includes a Sequence Listing.

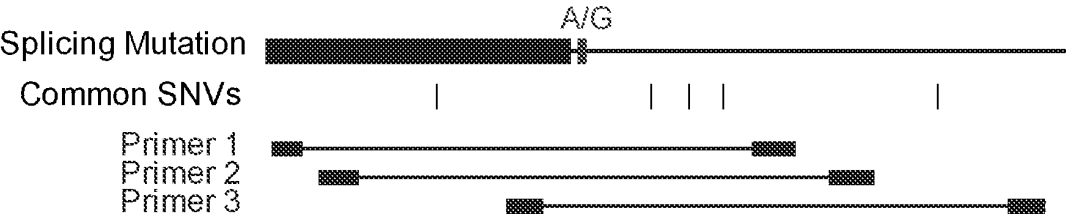


FIG. 1A

Adapter 1	Forward	Product 1	Reverse	Index 1	Adapter 2
Adapter 1	Forward	Product 2	Reverse	Index 2	Adapter 2
Adapter 1	Forward	Product 3	Reverse	Index 3	Adapter 2

FIG. 1B

Adapter 1	Forward	Product 1	Reverse	Index 1	UMI	Adapter 2
Adapter 1	Forward	Product 2	Reverse	Index 2	UMI	Adapter 2
Adapter 1	Forward	Product 3	Reverse	Index 3	UMI	Adapter 2

FIG. 1C

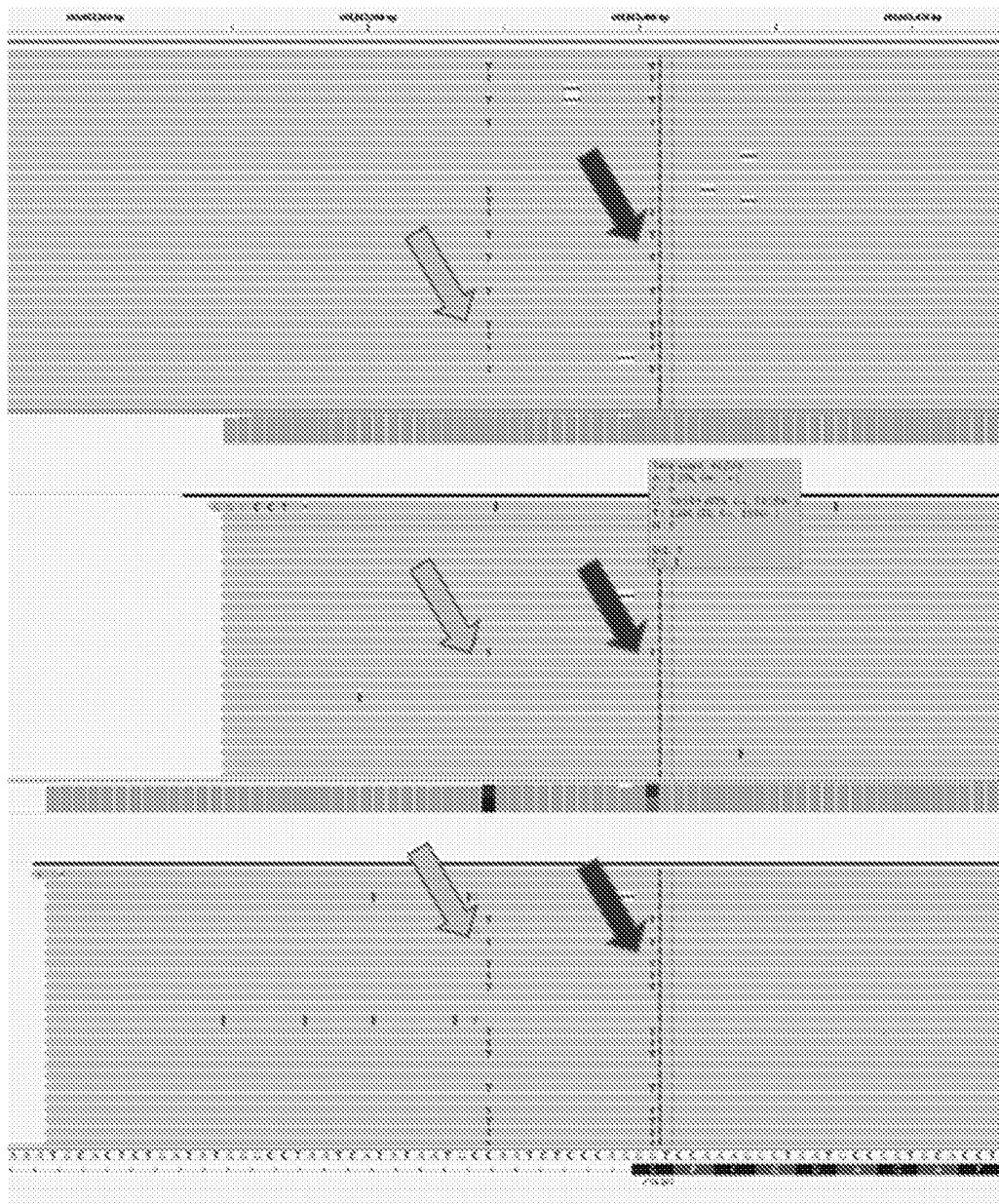


FIG. 2

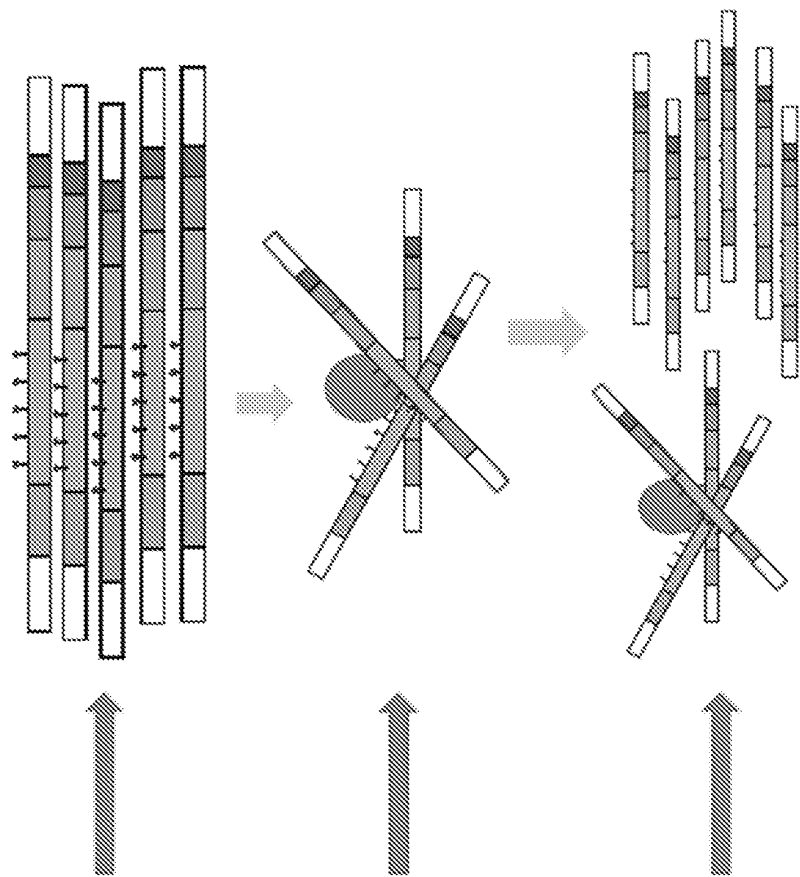


FIG. 3

- Low Cycle PCR
- Biotin labeling ~ 1/3 of all dCTPs
 - Independent reactions per primer
 - 25 to 50ng of input DNA
 - < 1ng of cDNA
- Purification of Sample
- Capture Amplified Sample on Streptavidin Beads
 - Wash to remove previous PCR mixture
- On-bead enrichment
- Custom primers
 - New product in solution
 - Maintains UMI and sample indexing

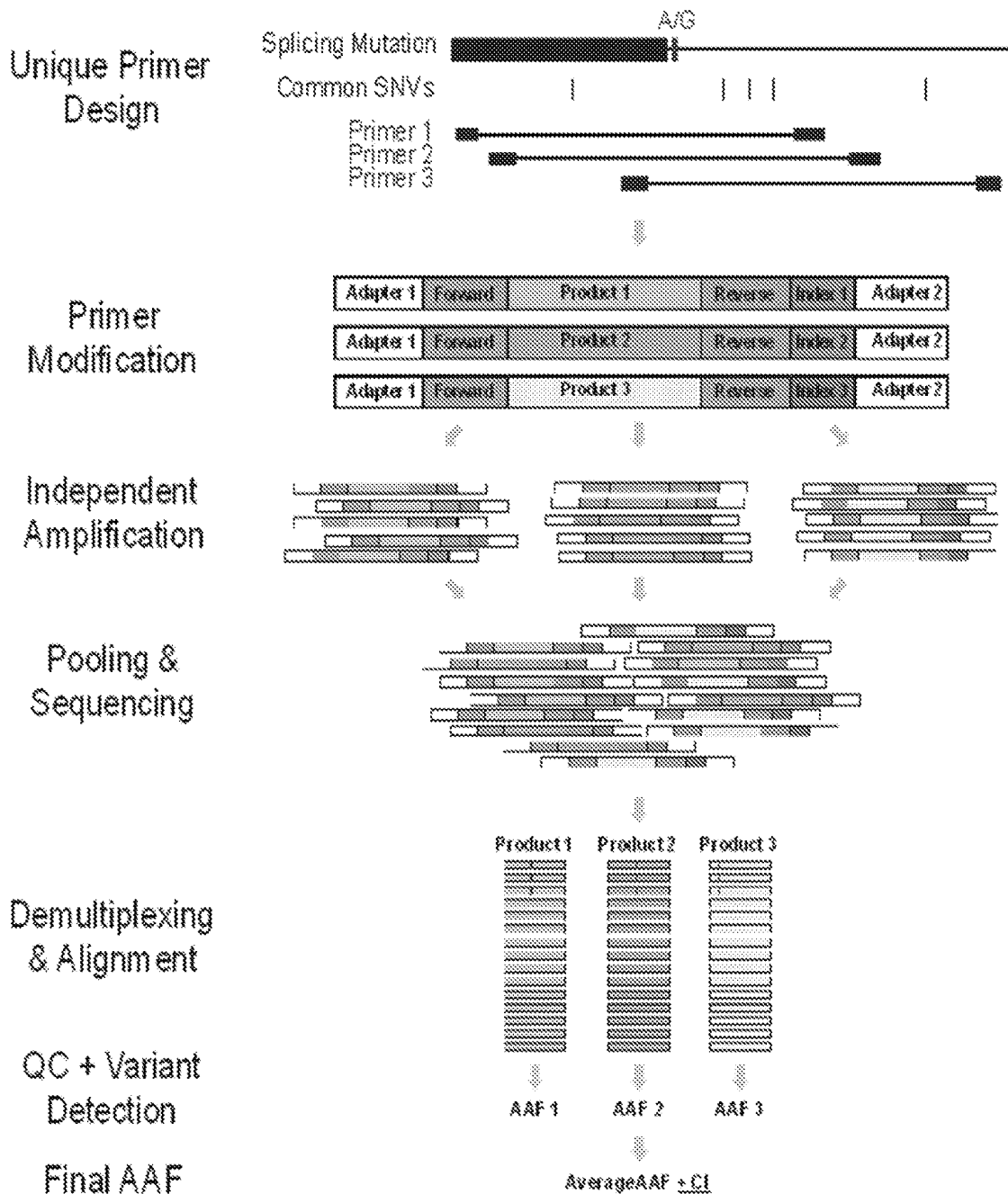
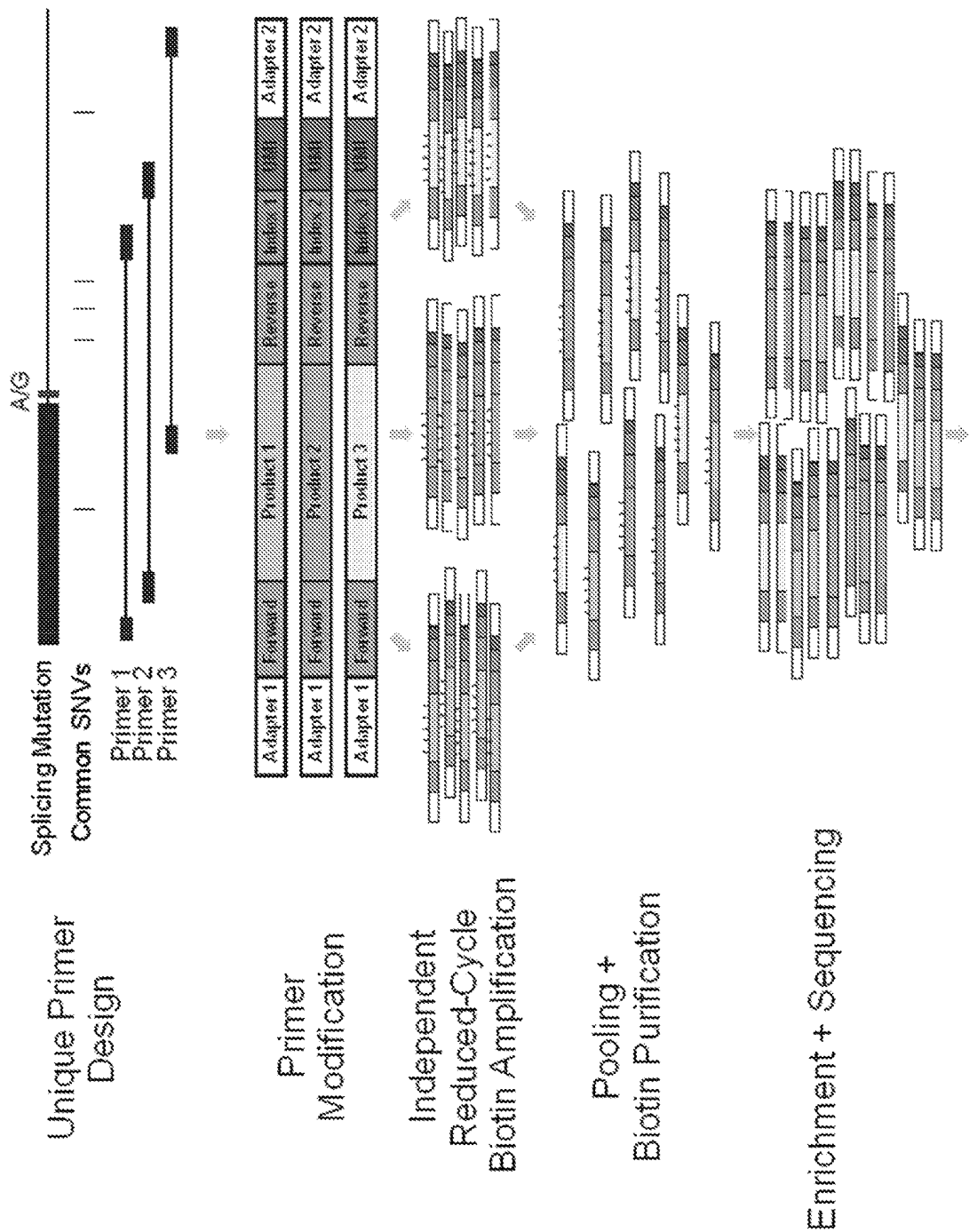


FIG. 4



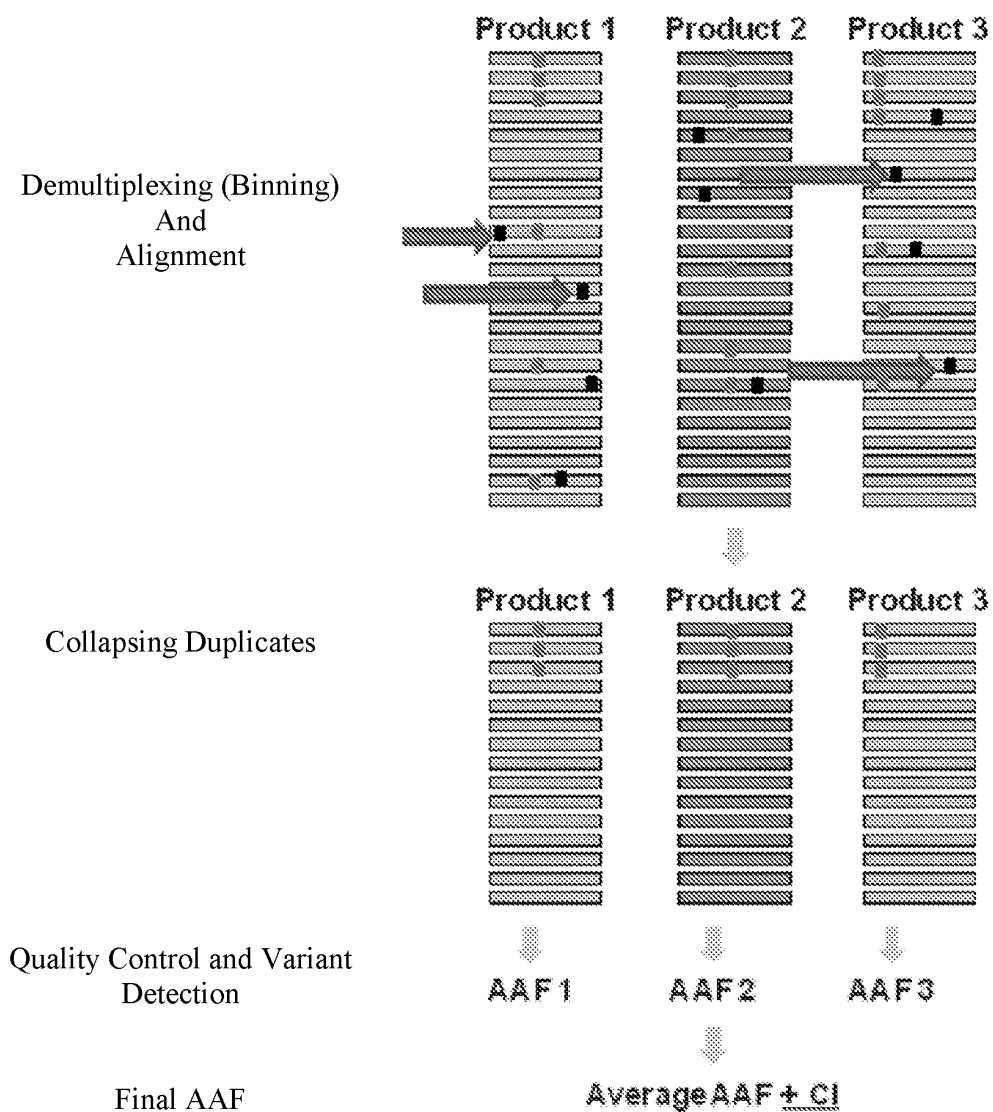


FIG. 5A (Continued)

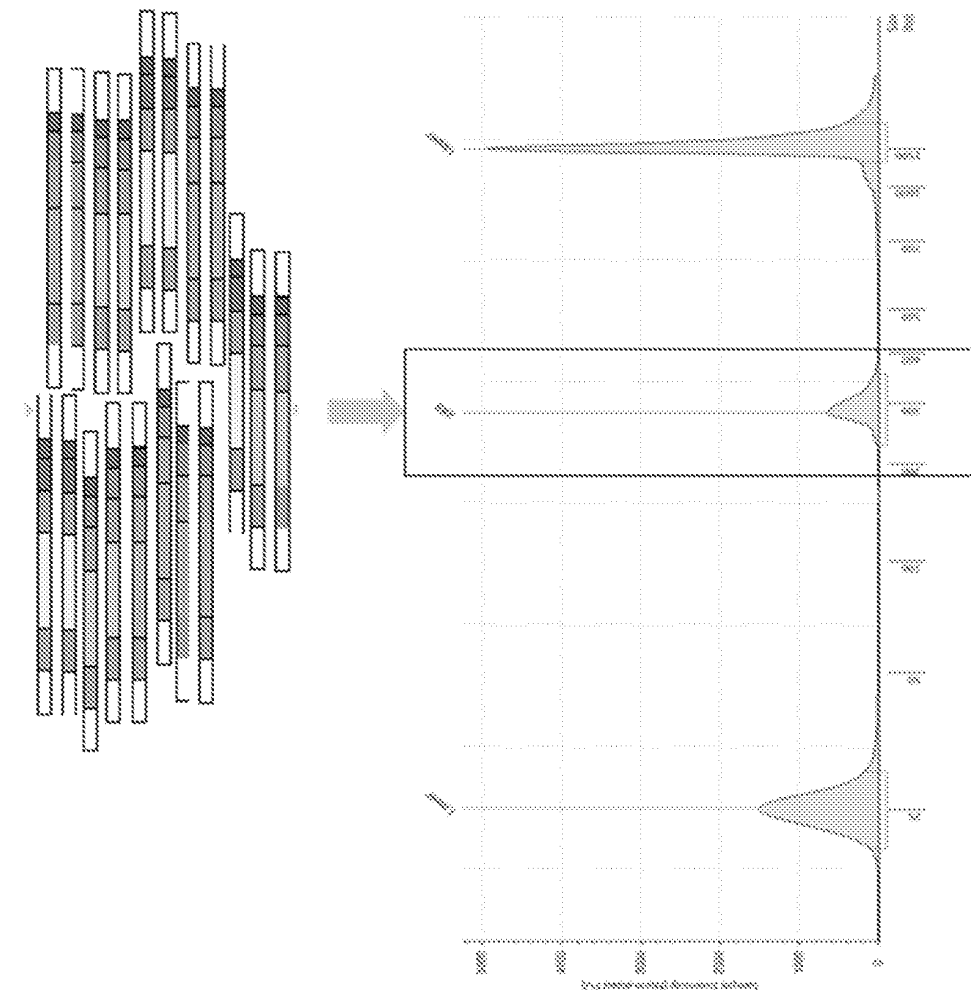


FIG. 5B

Purification:

- Remove biotin capture sample
- Pool all samples
- Purification using standard methods to remove leftover primers (MagJet or AmpureXP)

Bioanalyzer / TapeStation QC:

- Confirm all primer removed
- Single peak at appropriate size
- Concentration > 0.5ng/ul

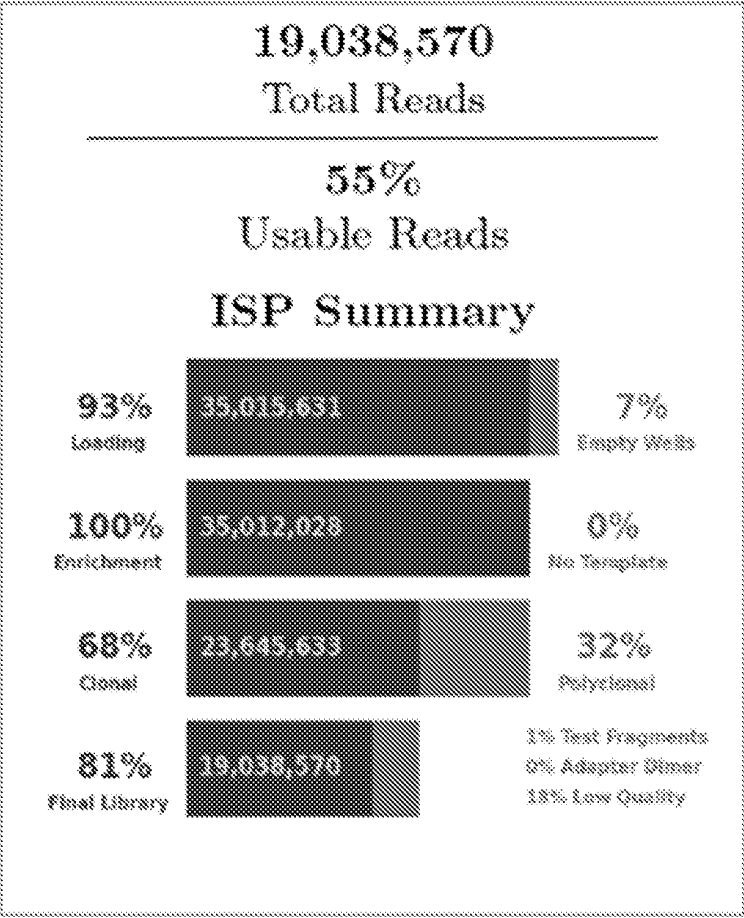


FIG. 6

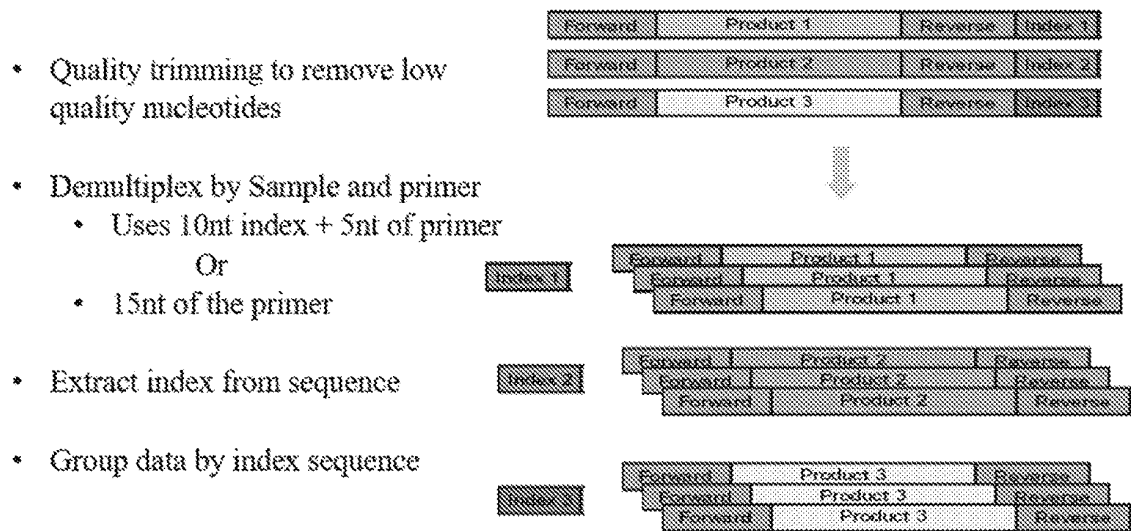


FIG. 7

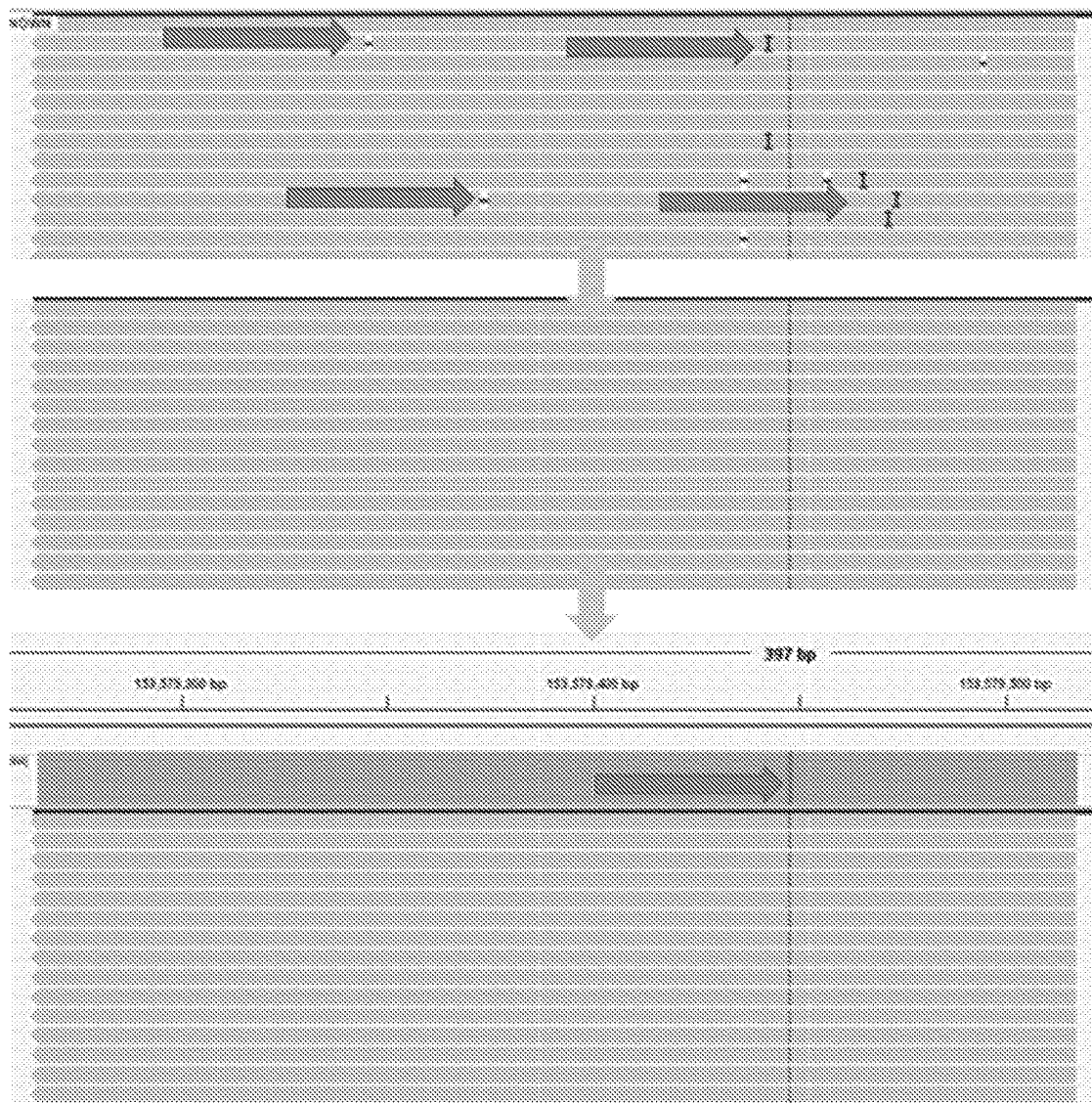


FIG. 8

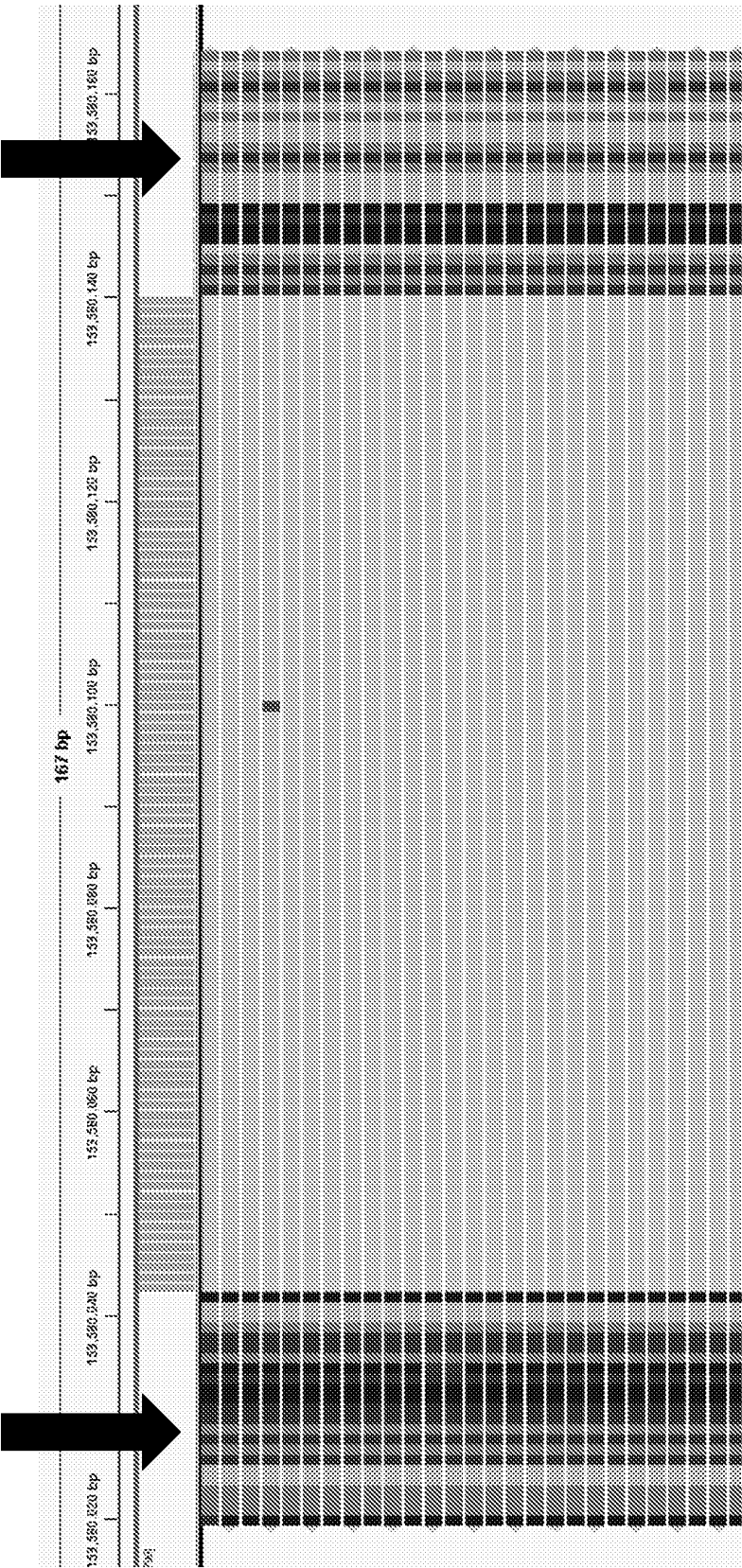


FIG. 9

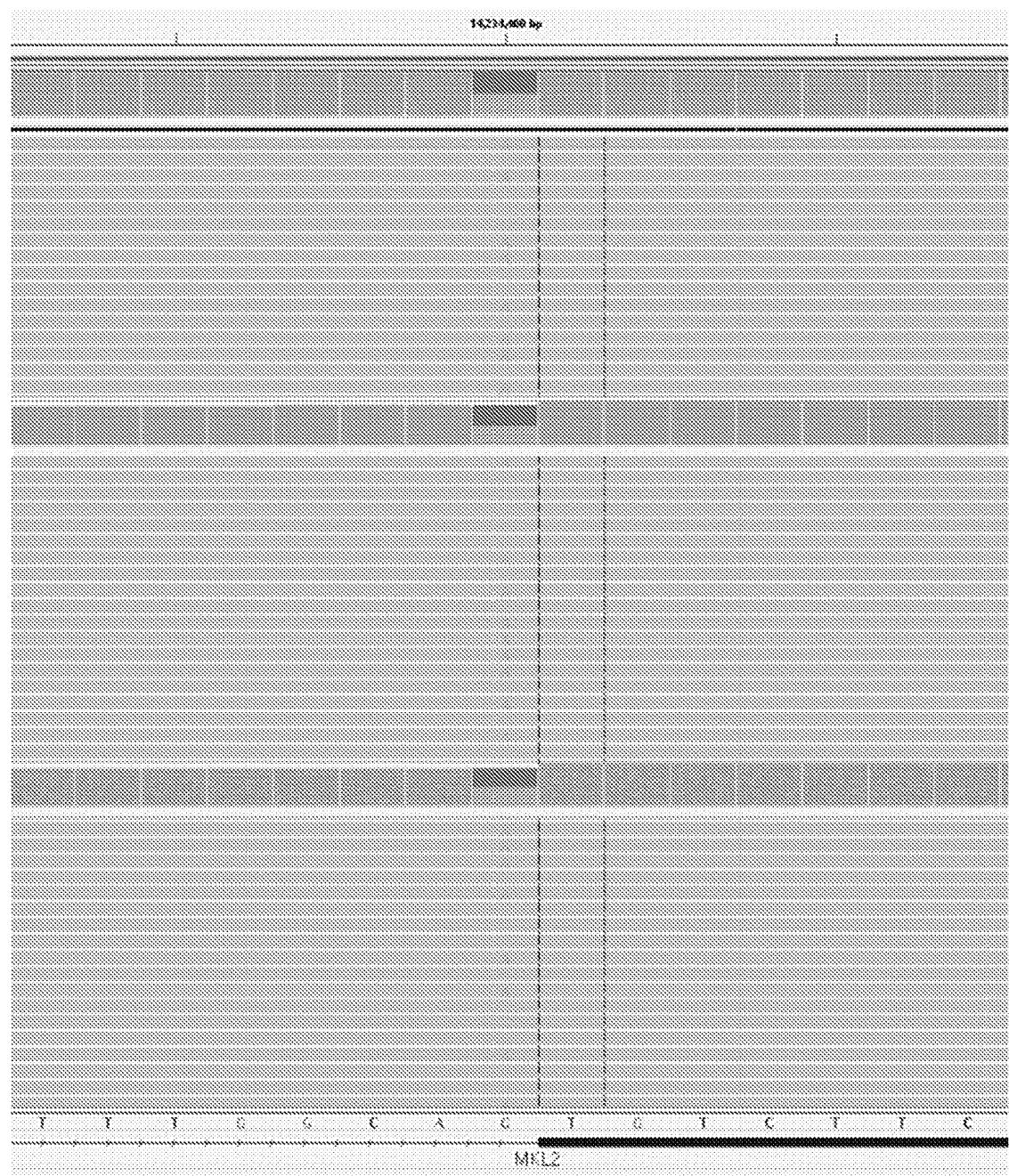


FIG. 10A

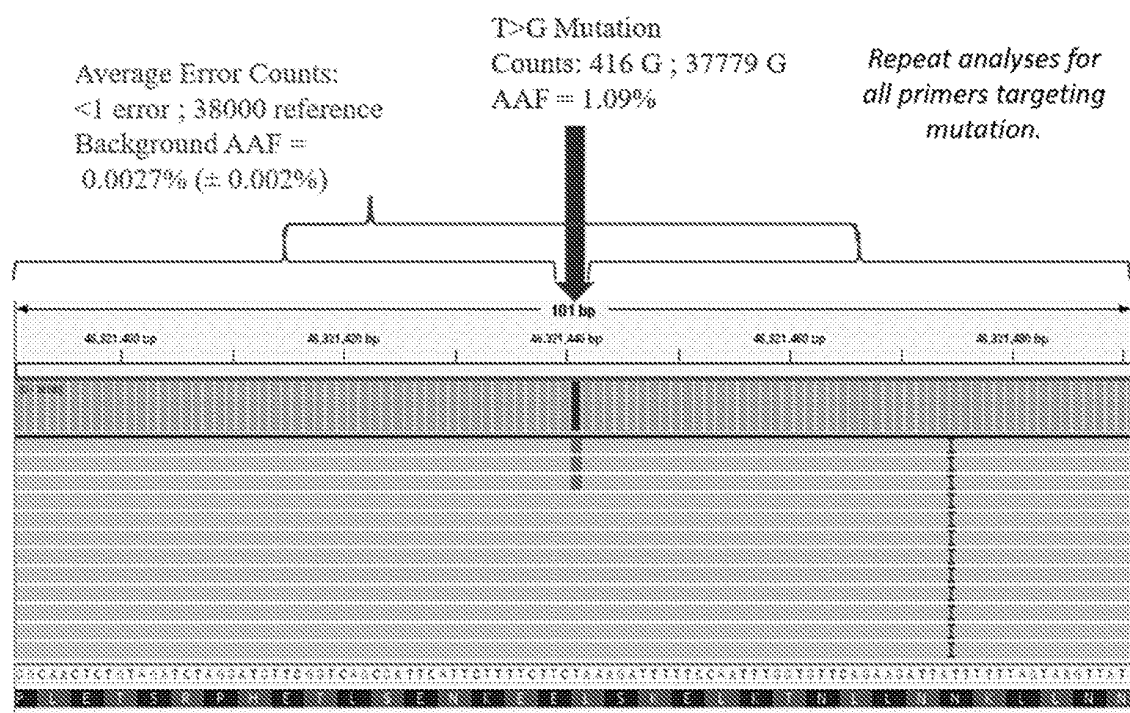


FIG. 10B

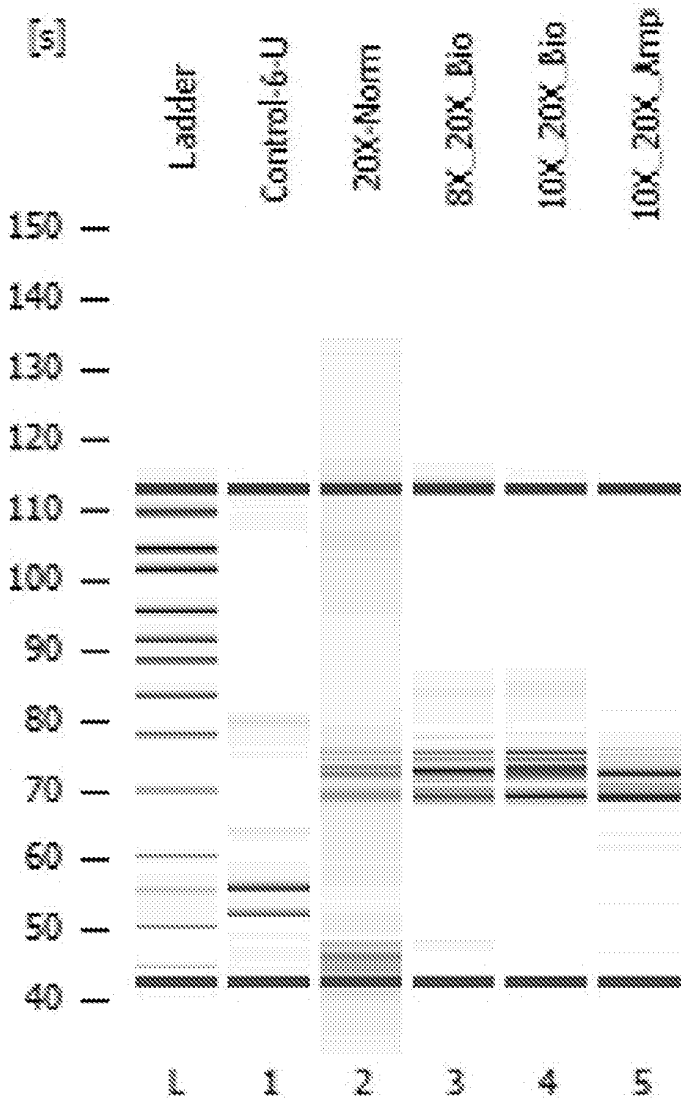


FIG. 11A

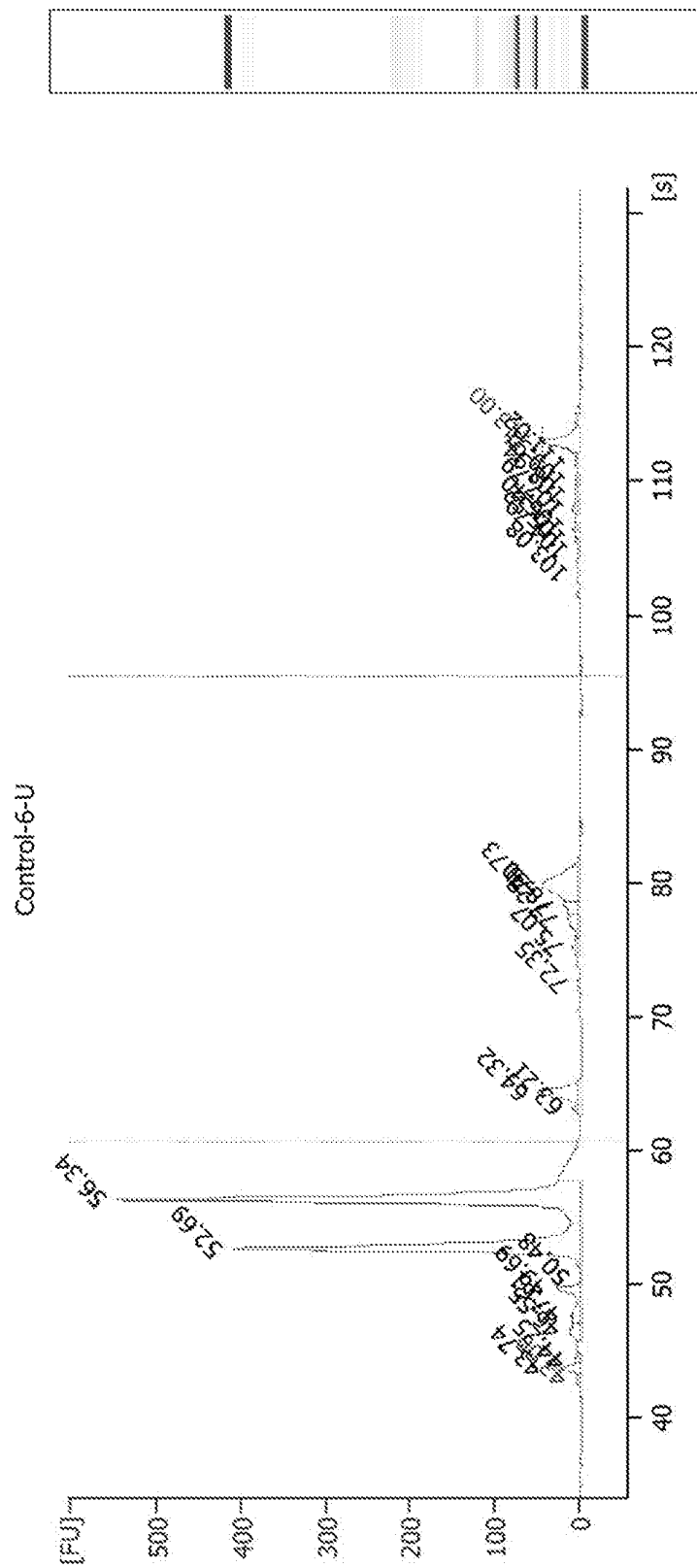


FIG. 11B

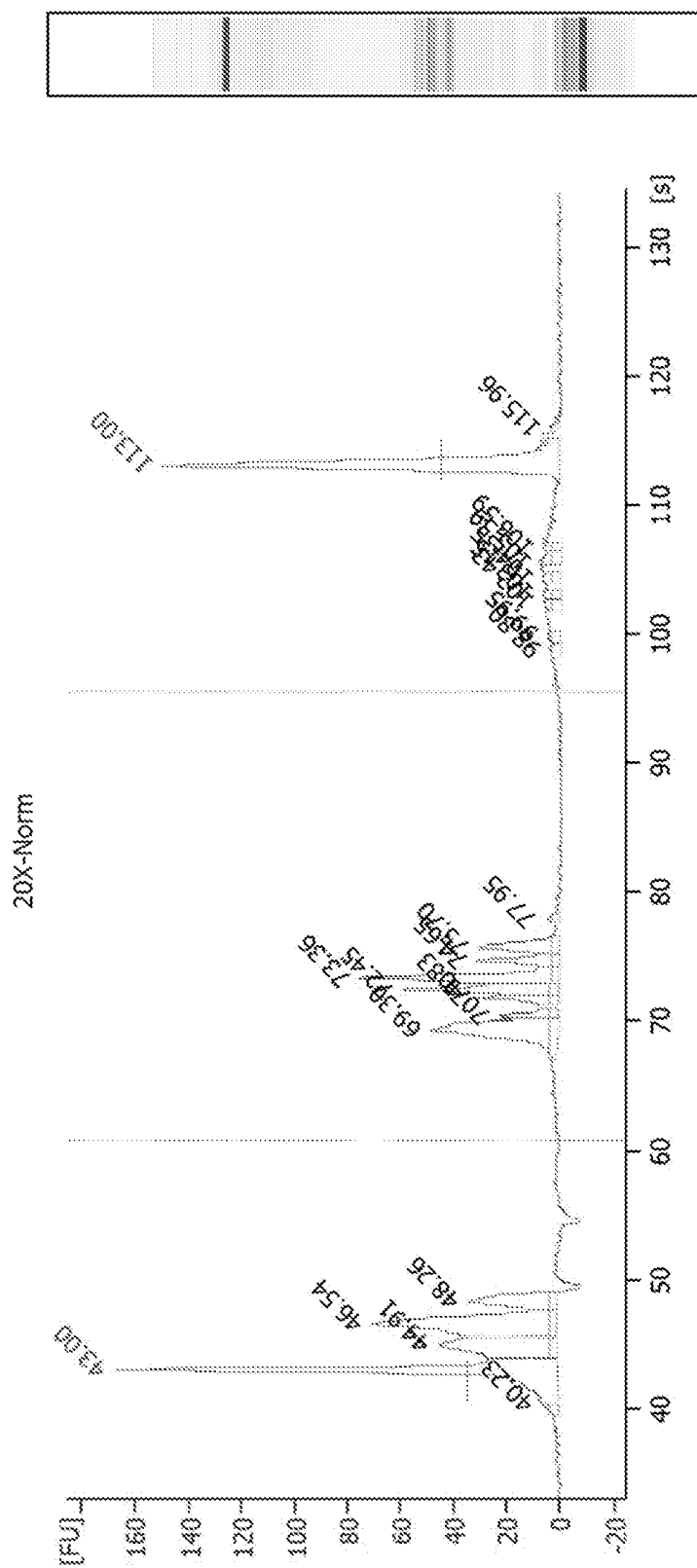


FIG. 11C

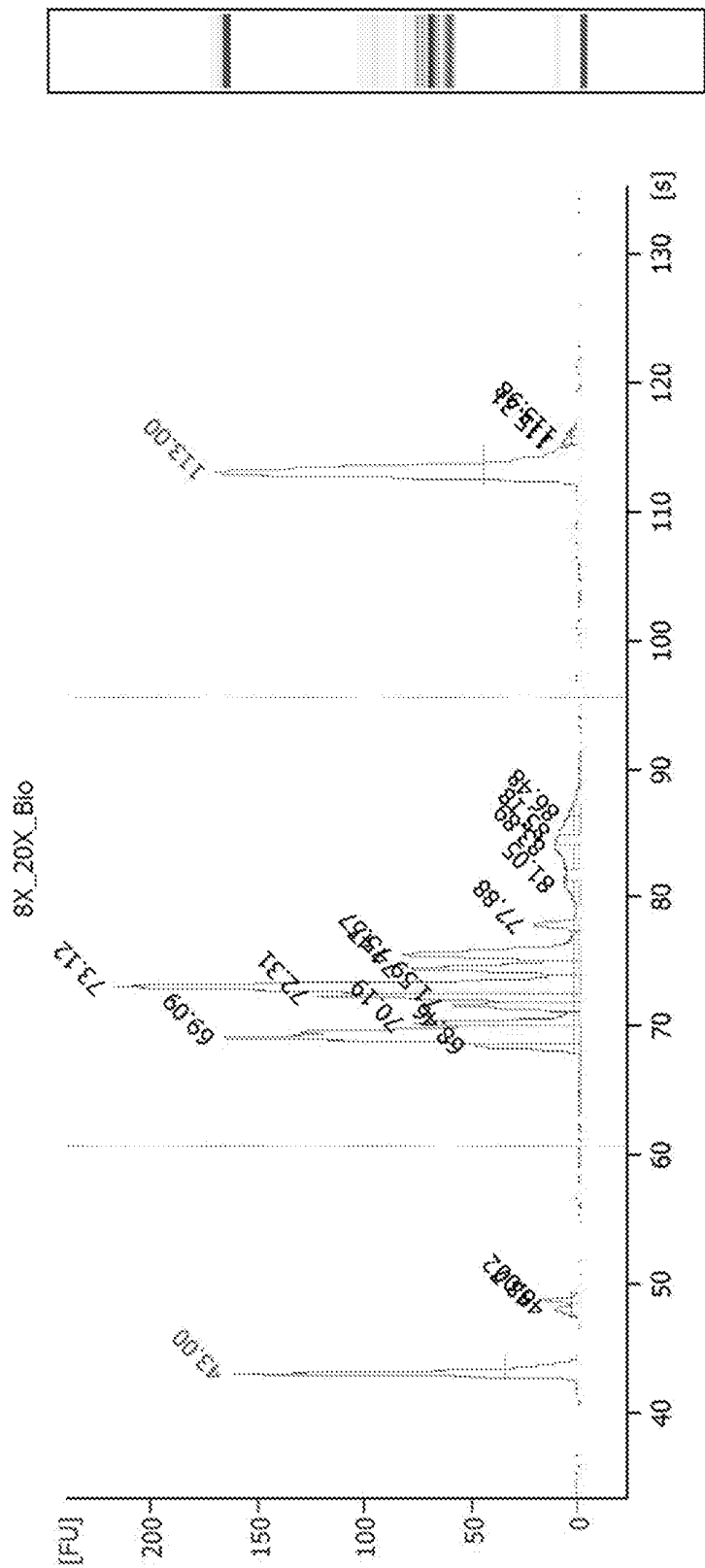


FIG. 11D

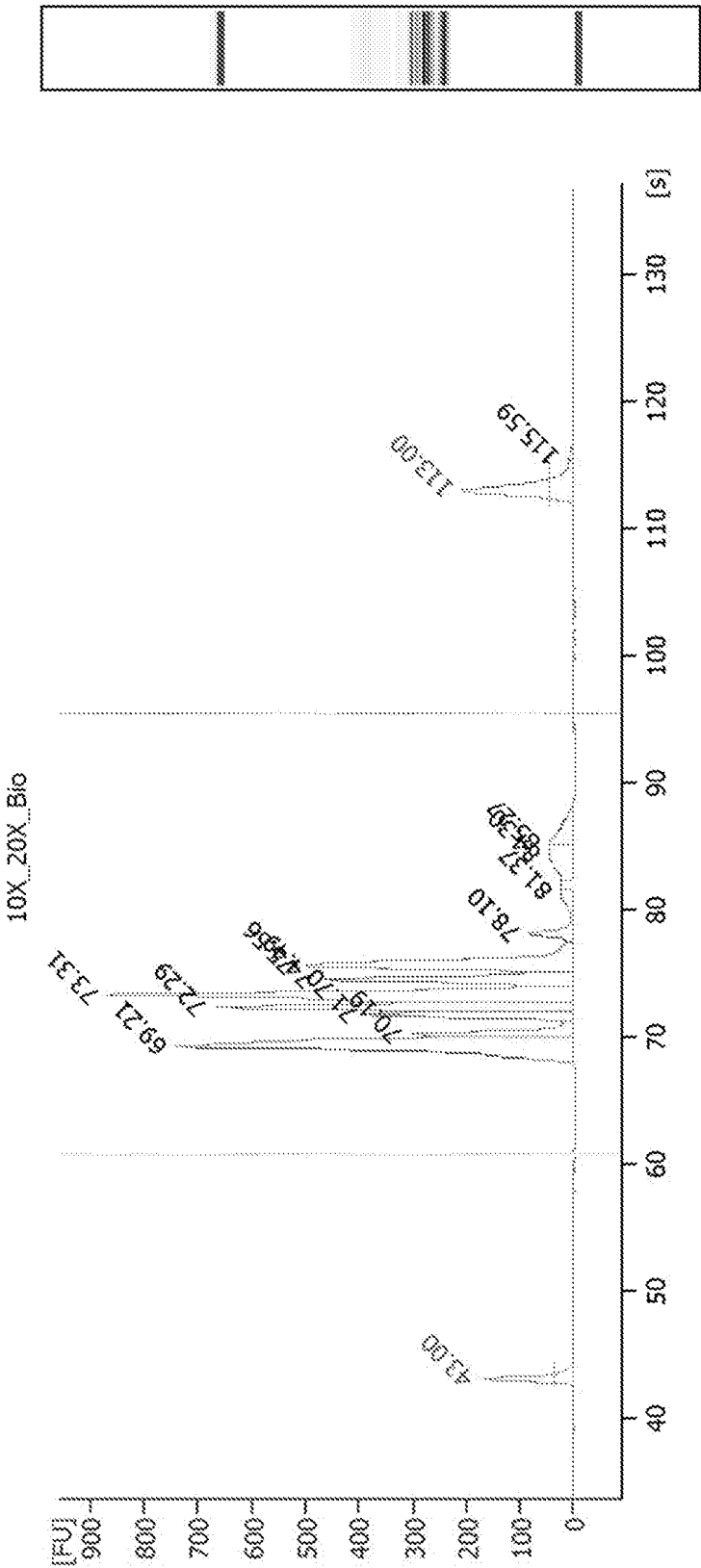


FIG. 11E

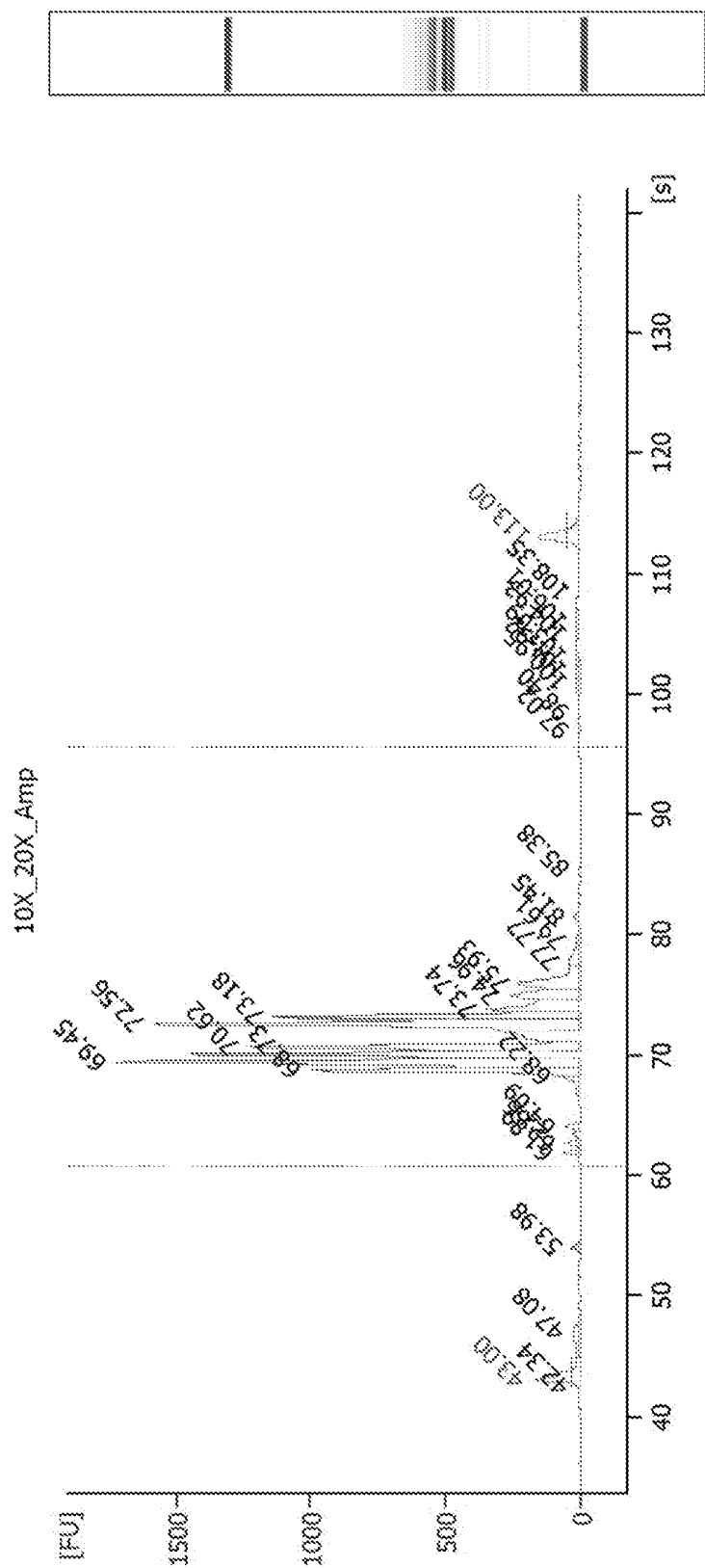


FIG. 11F

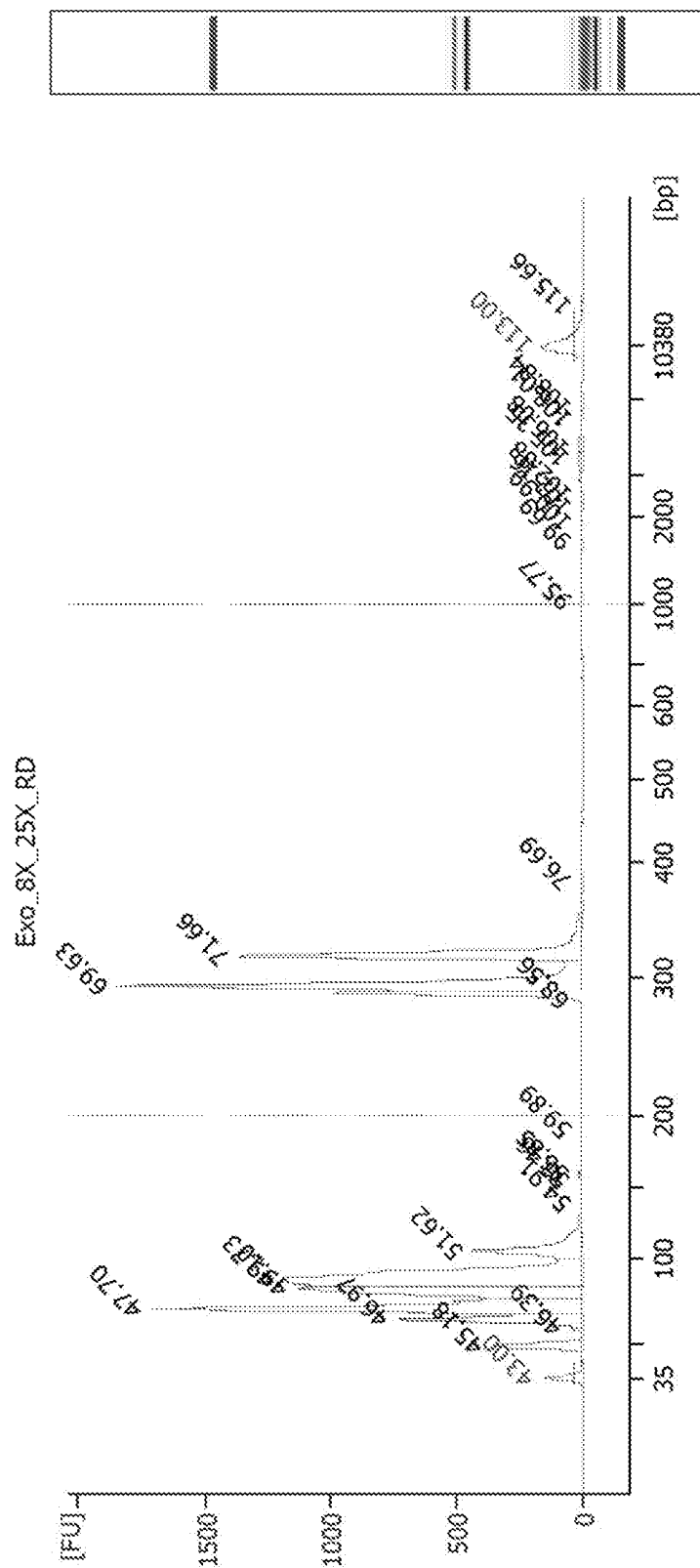


FIG. 11G

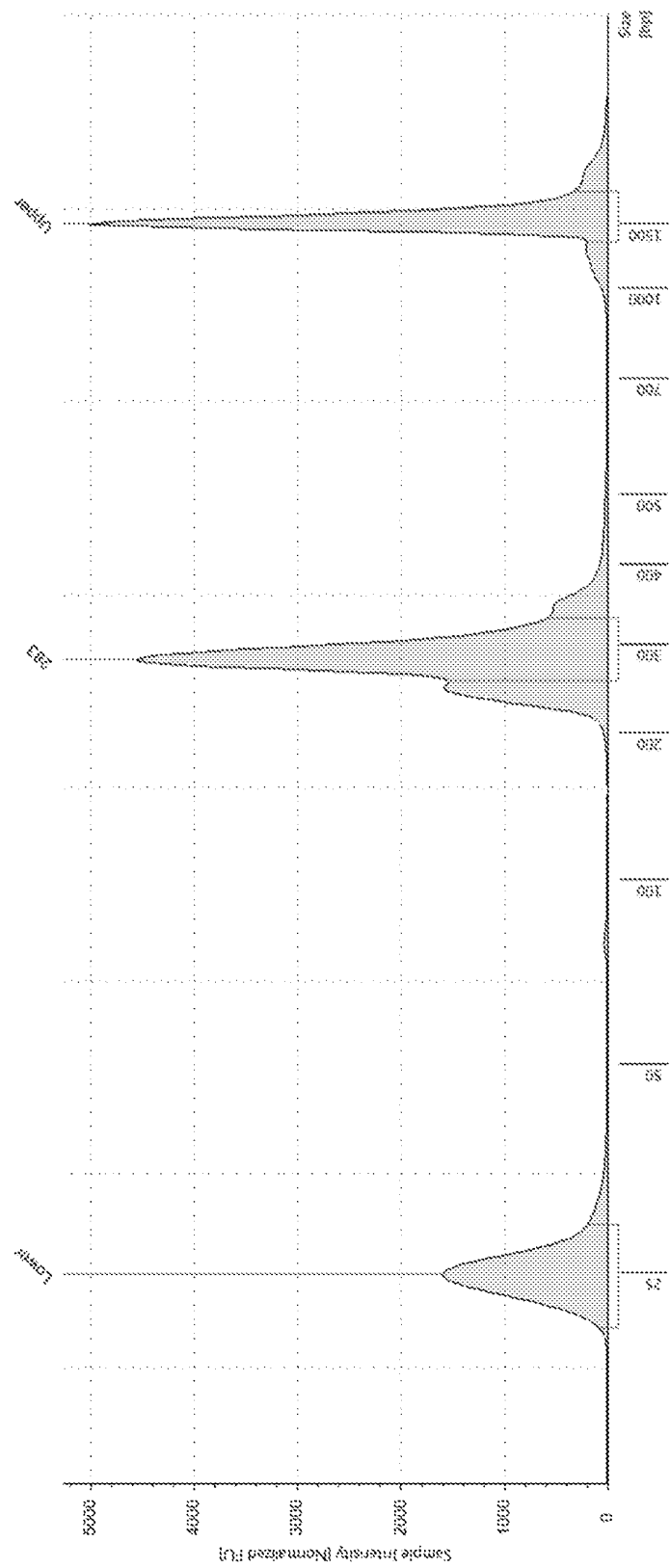


FIG. 12

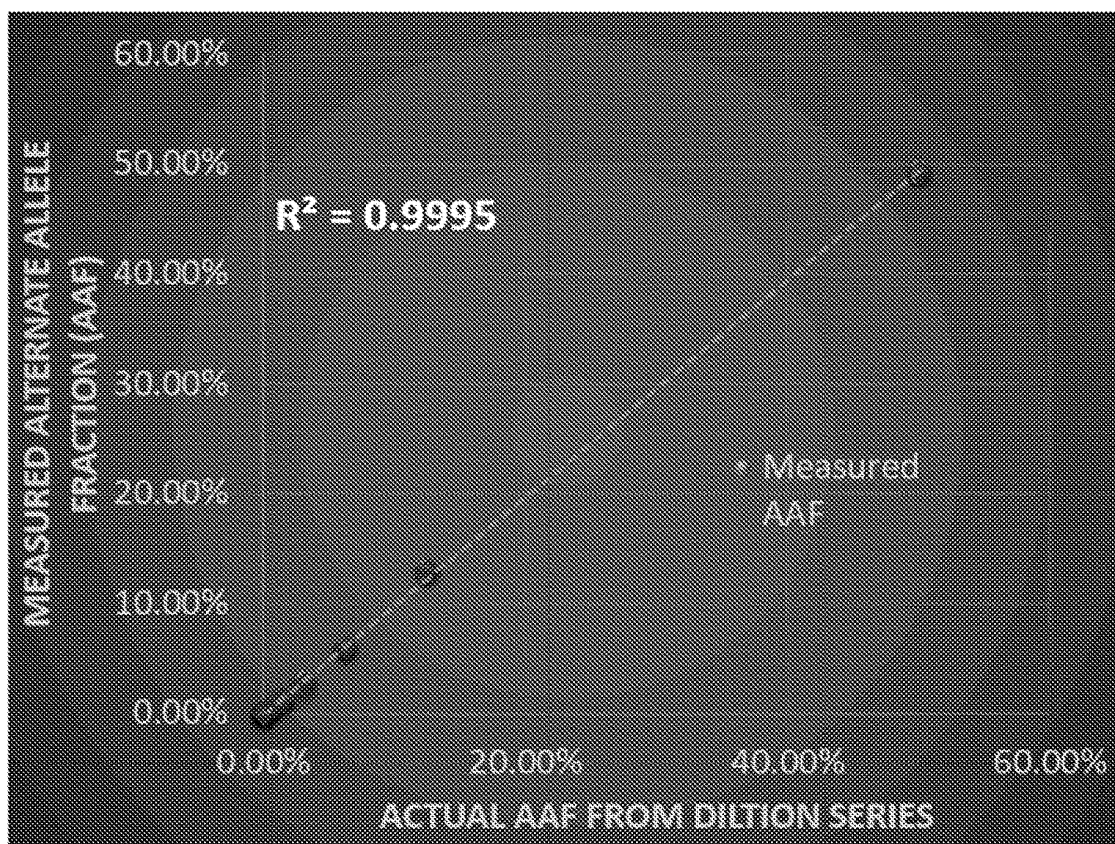


FIG. 13

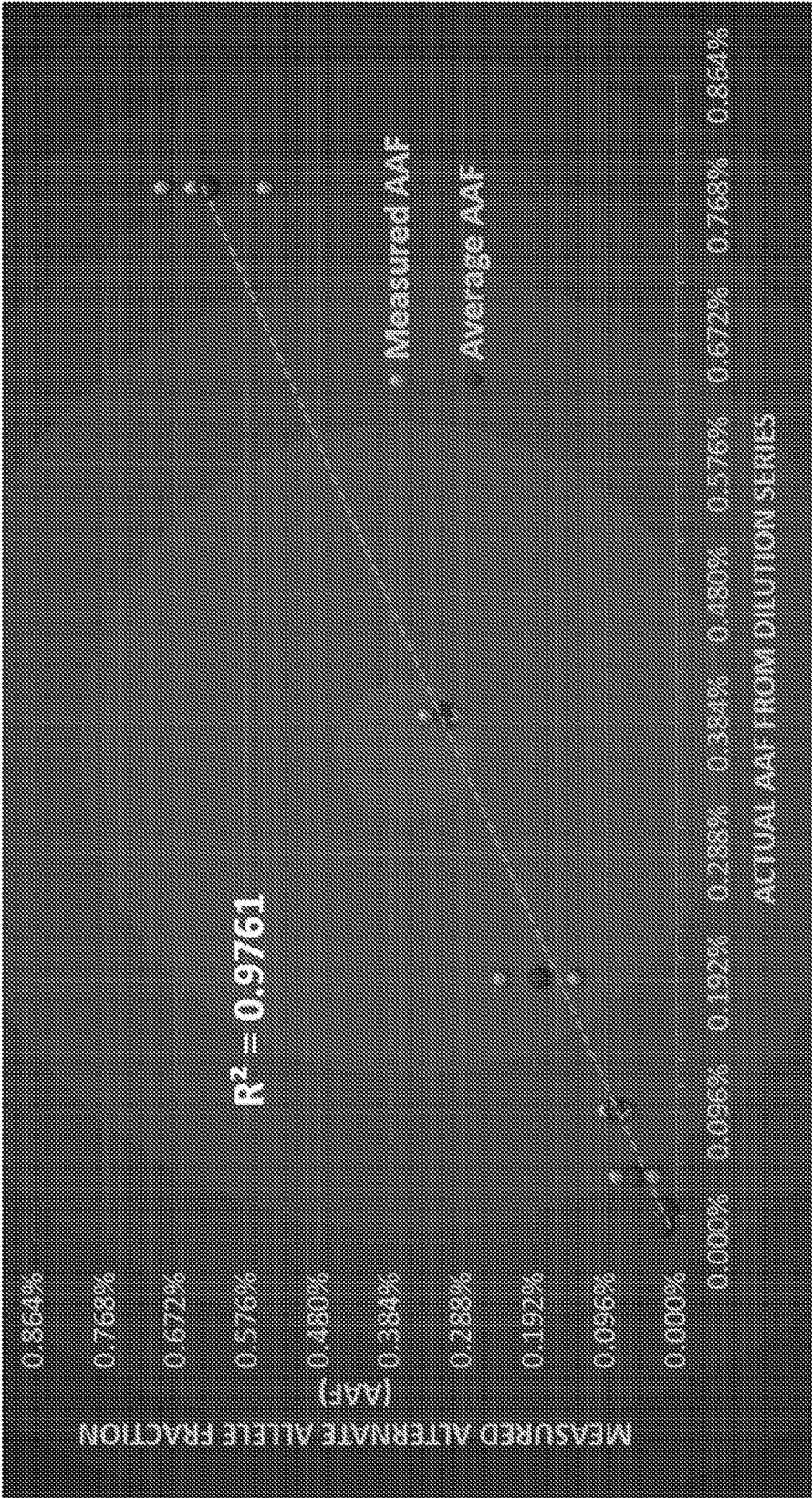


FIG. 14

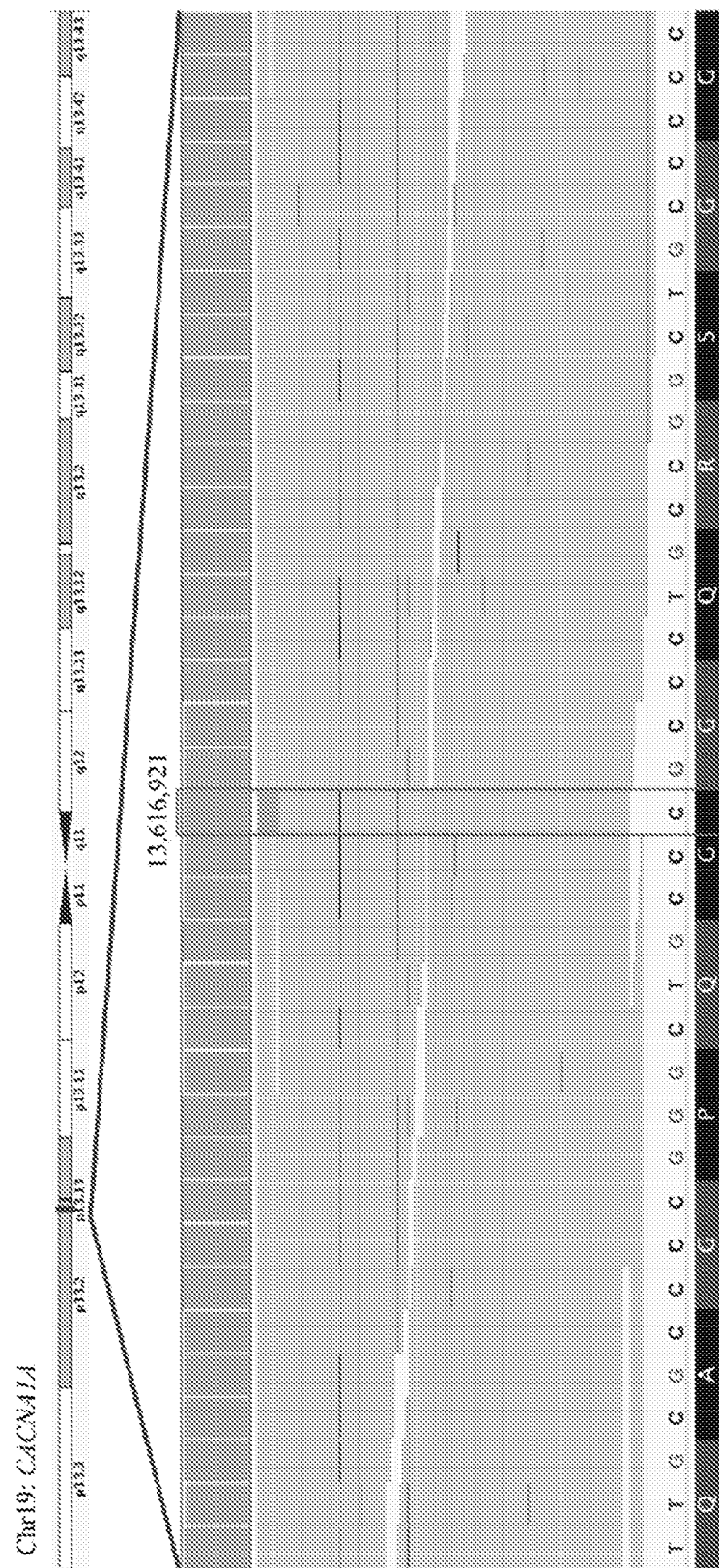


FIG. 15

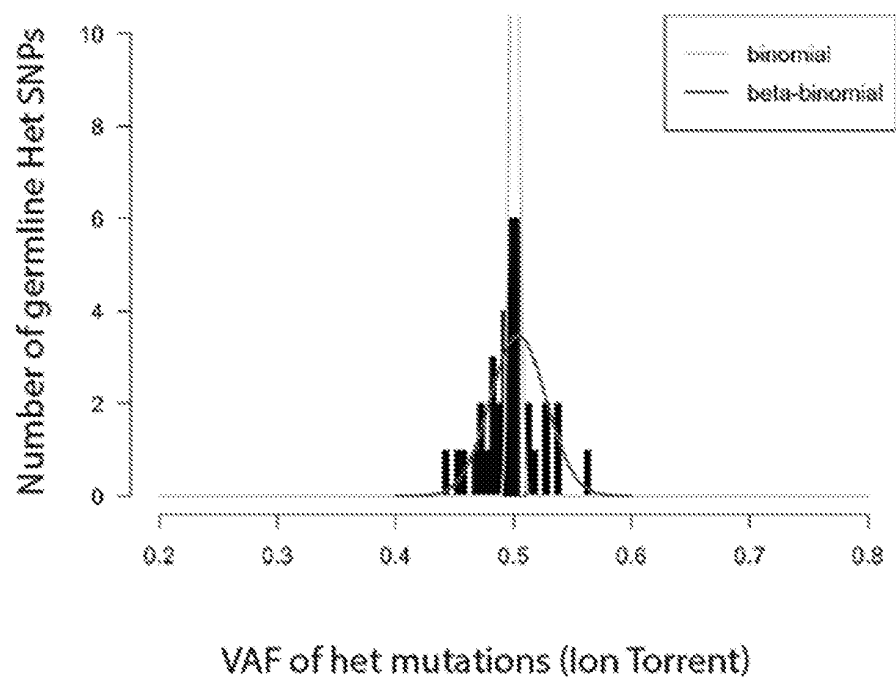


FIG. 16

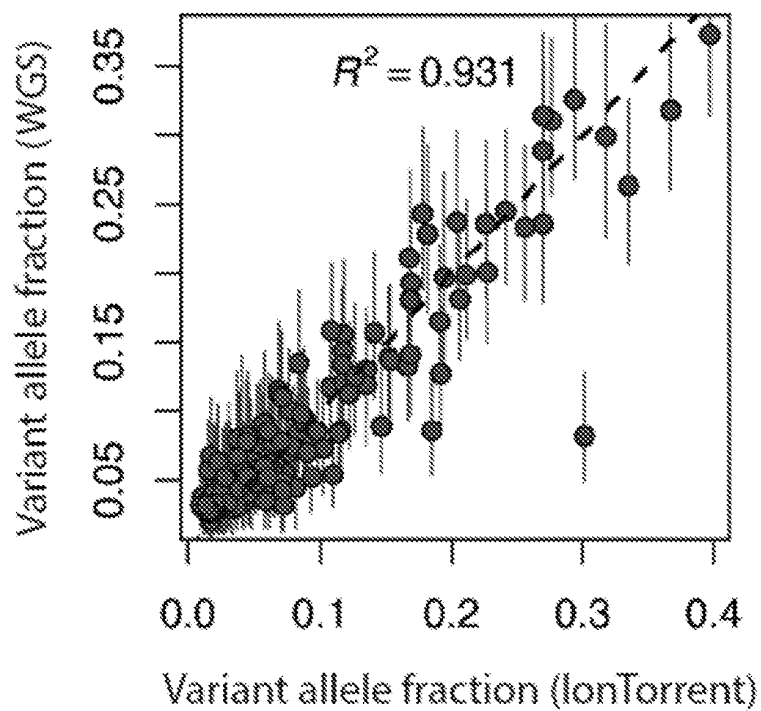


FIG. 17A

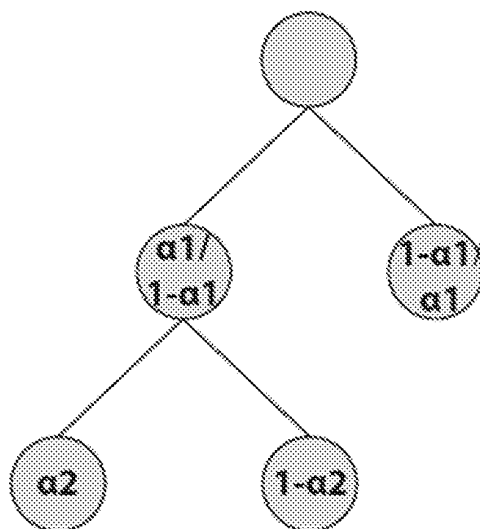


FIG. 17B

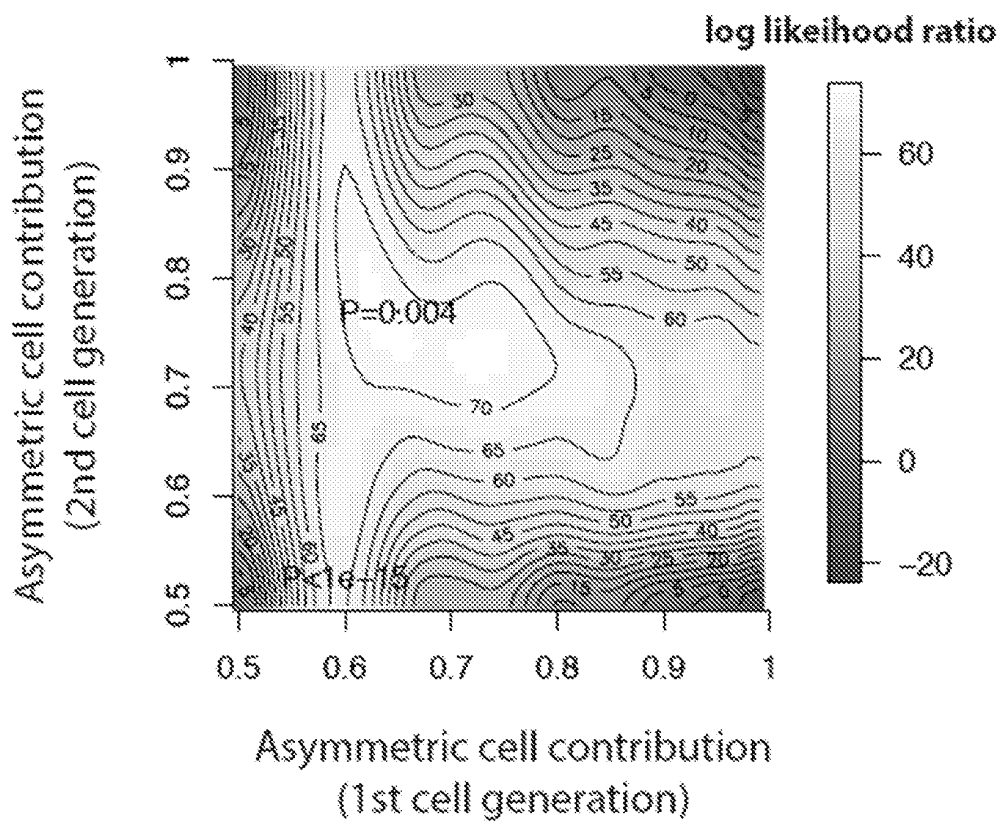


FIG. 17C

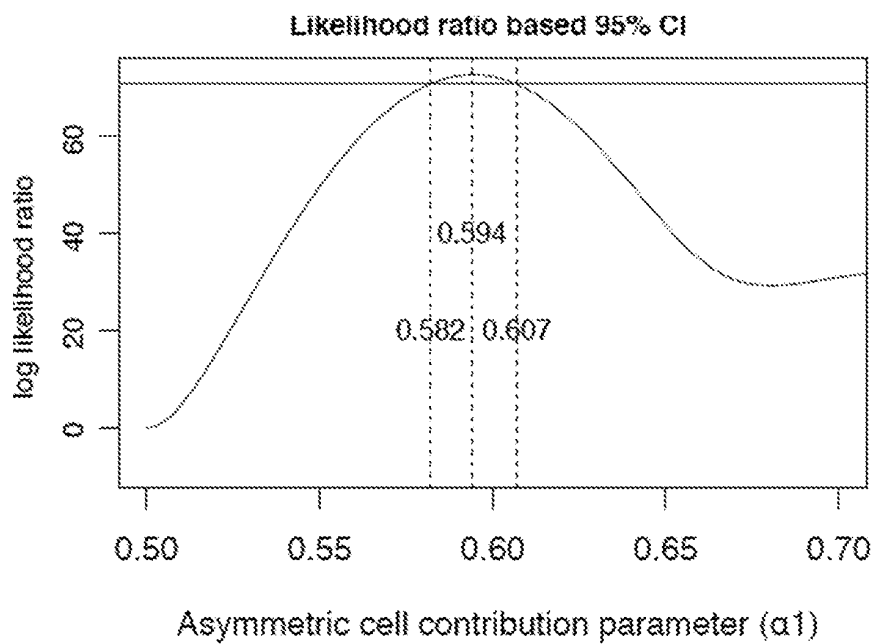


FIG. 17D

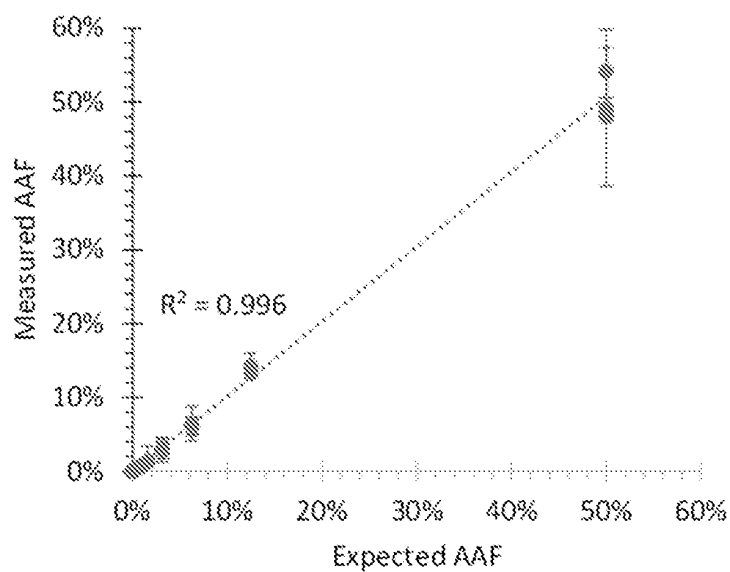


FIG. 18A

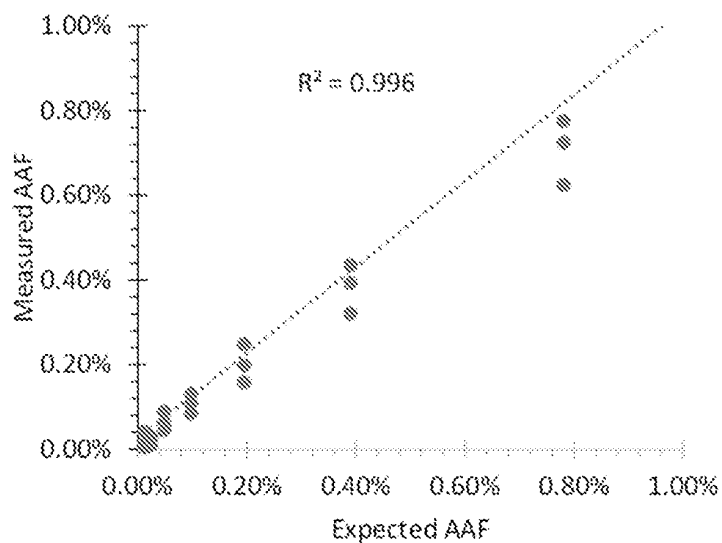


FIG. 18B

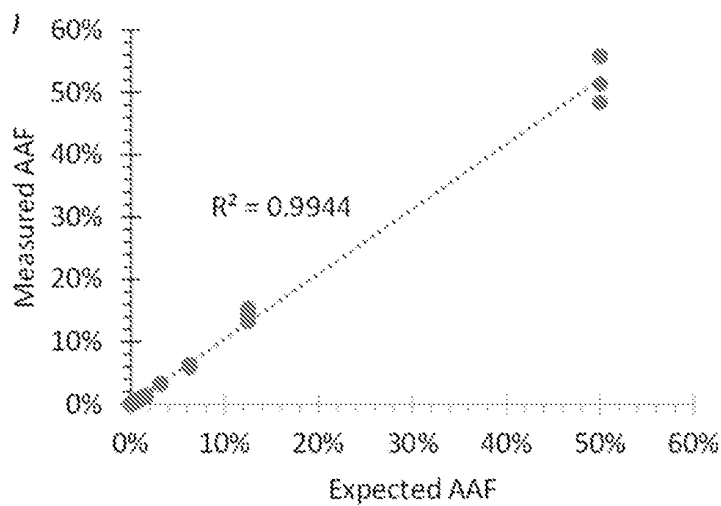


FIG. 18C

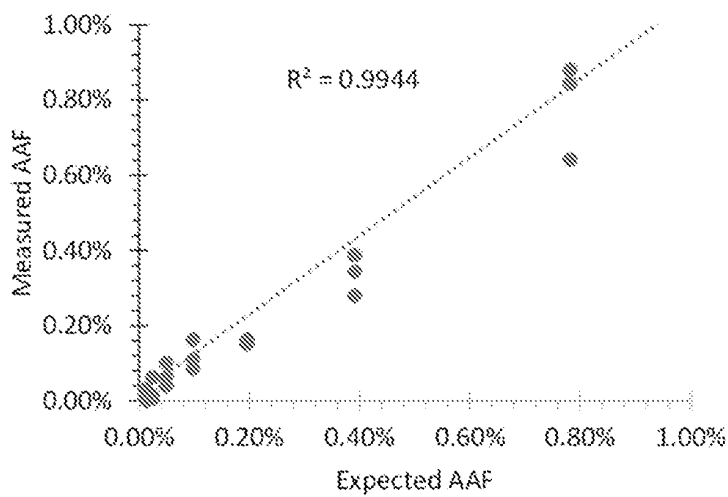


FIG. 18D

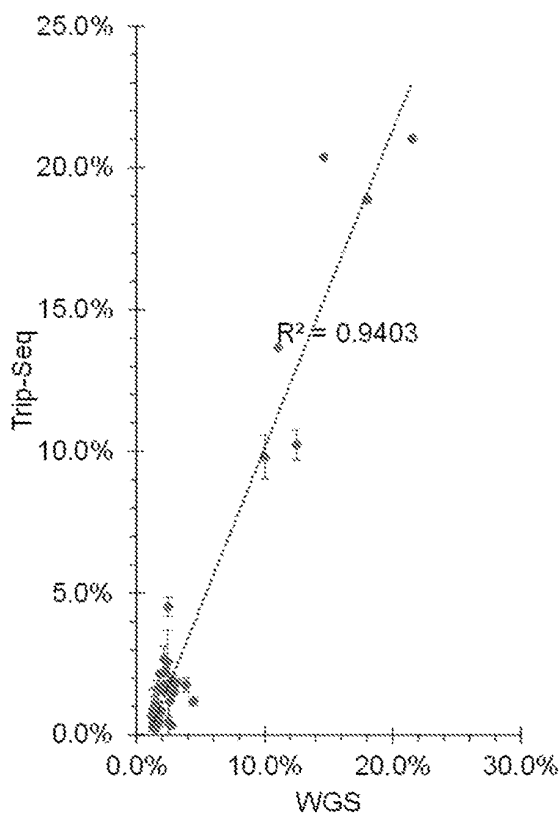


FIG. 19A

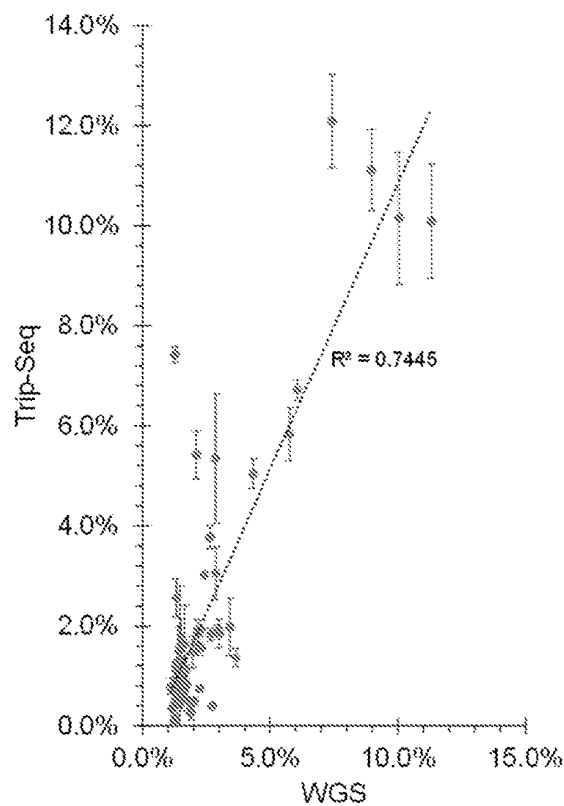


FIG. 19B

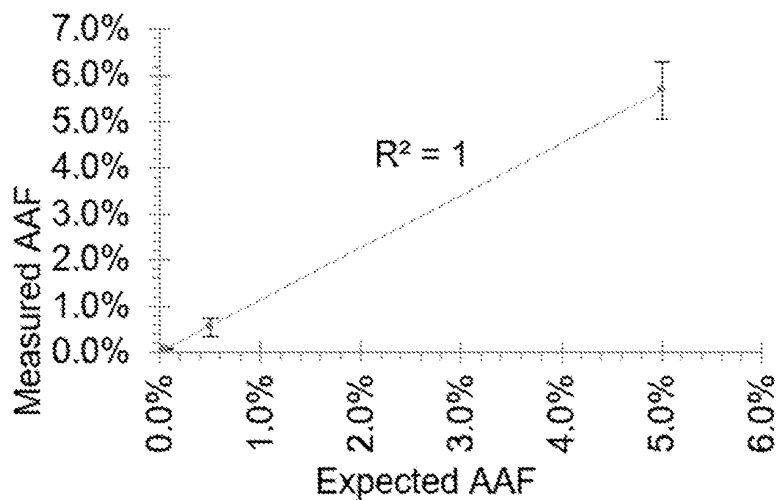


FIG. 19C

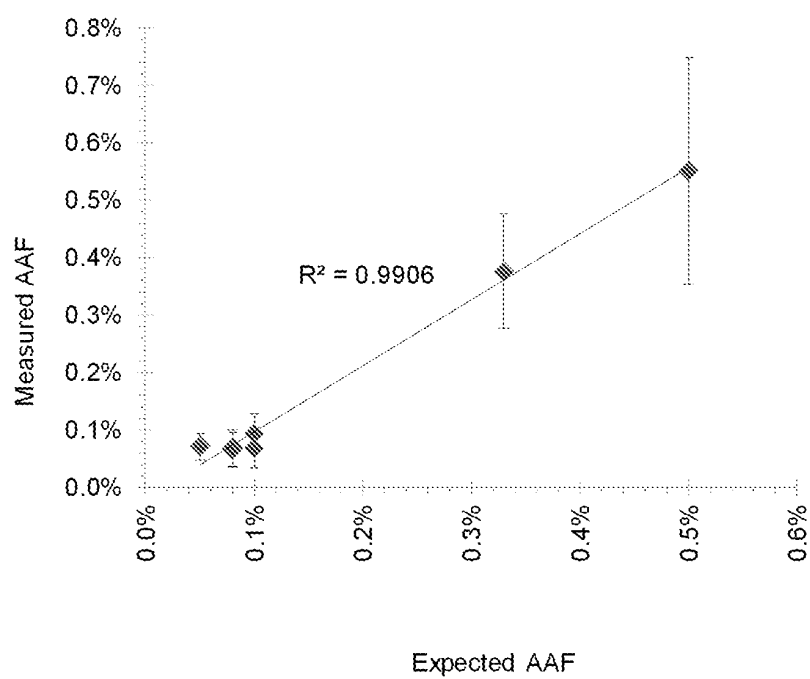


FIG. 19D

COST-EFFECTIVE DETECTION OF LOW FREQUENCY GENETIC VARIATION

CROSS-REFERENCE TO RELATED APPLICATION

[0001] This application claims the benefit of the following U.S. Provisional Application No. 62/799,671, filed Jan. 31, 2019, the entire contents of which are incorporated herein by reference.

STATEMENT OF RIGHTS TO INVENTIONS MADE UNDER FEDERALLY SPONSORED RESEARCH

[0002] This invention was made with government support under Grant Nos. R01NS032457 and U01MH106883 awarded by the National Institutes of Health. The government has certain rights in the invention.

BACKGROUND OF THE INVENTION

[0003] Traditional genetic sequencing methodologies, such as whole genome (WGS) and whole exome (WES), have focused on the important contribution of germline mutations that are present in all cells throughout the human body. However, recent studies have shown numerous examples of mutations occurring after fertilization (i.e. postzygotic mutations), which are only present in a fraction of the cells. Postzygotic mutations, or somatic mutations, have been heavily studied in cancers where clinical diagnostic testing for somatic mutations in tumor and blood samples are becoming a standard practice due to improved detection sensitivities when most cells in the sample carry a given mutation.

[0004] Beyond technical errors, an important consideration for skewed alternate allelic fraction (AAFs), false negatives, and false positives are allelic imbalances caused by inherent differences in the genome content around a mutation. These issues, such as additional mutations, repeat content, methylation, or copy number changes, can have dramatic impacts on AAFs, resulting in the commonly recognized issue of allelic dropout. To avoid allelic dropout, many methods avoid placing primers in areas with known genetic variation in the general population. However, these methods remain susceptible to allelic skewing from ultrarare or private alleles and other loci specific causes of allelic imbalance. Cost-effective methods are needed for the detection and characterization of rare alleles and other genetic variants.

SUMMARY OF THE INVENTION

[0005] As described below, the present disclosure features methods for detecting and quantifying genetic variants in a sample.

[0006] In one aspect of the present disclosure, a method is provided for determining alternate allele frequency, the method involves performing two or more parallel amplification reactions on a single sample, thereby generating overlapping amplicons, where each amplification reaction includes a unique pair of forward and reverse primers, where the forward or reverse primer includes an index sequence, and where the forward and reverse primers include different adapter sequences. The method also involves sequencing the overlapping amplicons to produce sequence reads, segregating the sequencing reads into bins by index sequence, and

detecting the presence or absence of one or more genetic variants within sequencing reads within a bin, where the frequency of detection of the variant determines the alternate allele frequency.

[0007] Another aspect provides a method for determining alternate allele frequency, the method involves a) performing three amplification reactions on a single sample, thereby generating three overlapping amplicons, where each amplification reaction includes a unique pair of forward and reverse primers, where each primer includes a nucleic acid sequence complementary to a portion of a target nucleic acid sequence, where the forward or reverse primer includes an index sequence, where the forward and reverse primers include different adapter sequences at or near the 5' terminus of the primer and upstream of the sequence complementary to the target, and where at least one adapter sequence is complementary to a nucleic acid sequence used in sequencing; b) sequencing the overlapping amplicons to produce sequence reads; c) segregating the sequencing reads into bins by index sequence; and d) detecting the presence or absence of one or more genetic variants within sequencing reads within a bin, where the frequency of detection of the variant determines the alternate allele frequency.

[0008] Another aspect of the present invention provides a method for method for determining alternate allele frequency, the method involving a) performing three amplification reactions on a single sample, thereby generating three overlapping amplicons, where each amplification reaction includes a unique pair of forward and reverse primers, where the forward or reverse primer comprises an index sequence and/or a unique molecular identifier (UMI); and each primer includes i. a nucleotide sequence complementary to a portion of a target nucleic acid sequence; ii. an adapter at or near its 5' terminus, where the adapter is upstream of the sequence complementary to the target and wherein the forward and reverse primers include different adapter sequences, and where at least one adapter sequence is complementary to a nucleic acid sequence used in sequencing; b) sequencing the overlapping amplicons to produce sequence reads; c) segregating the sequencing reads into bins by index sequence; d) detecting the UMI and removing duplicate reads from the bin, where the detecting can be simultaneous with step c or subsequent to step c; and e) detecting the presence or absence of one or more genetic variants within sequencing reads within a bin, where the frequency of detection of the variant determines the alternate allele frequency.

[0009] In some embodiments, the methods disclosed herein further involve pooling the amplicons prior to sequencing. In some embodiments of the methods disclosed herein, sequencing the amplicons involves contacting the amplicons with a nucleic acid complementary to the adapter sequence. In some embodiments, the amplicons include a nucleotide having a label, and in some embodiments, the label is biotin. In some embodiments, the methods disclosed herein also involve contacting the label with a capture agent that specifically binds the label. In some embodiments, the methods also involve enzymatically digesting the primers. In some embodiments of the present disclosure, the methods also involve amplifying the amplicons, thereby generating enriched populations of amplicons. In some embodiments, the genetic variation to be detected is known or unknown. In some embodiments, the genetic variant has an alternate allele fraction of at least 0.1%. In some embodiments, the

genetic variant has an alternate allele fraction of at least 0.025%. In some embodiments, the genetic variant is a mosaic variant. In some embodiments, detection of the genetic variant identifies the presence of a disease or a predisposition to a disease in a subject from whom the sample was derived. In some embodiments, the disease is cancer. In some embodiments, the sample includes circulating tumor cells or cell free DNA. In some embodiments, the genetic variant originated from a somatic event or a germline event. In some embodiments, the alternate allele frequency is compared to the allele frequency of a reference sample to determine if the subject's disease is progressing, regressing, or in remission. In some embodiments, the methods further involve averaging the alternate allele frequencies determined for each bin. In some embodiments, the methods further involve determining the error rate of the nucleic acid sequences flanking the alternate allele.

[0010] Methods defined by the present disclosure were performed in connection with the examples provided below. Other features and advantages of the disclosure will be apparent from the detailed description and from the claims.

Definitions

[0011] Unless defined otherwise, all technical and scientific terms used herein have the meaning commonly understood by a person skilled in the art to which this disclosure relates. The following references provide one of skill with a general definition of many of the terms used in this disclosure: Singleton et al., Dictionary of Microbiology and Molecular Biology (2nd ed. 1994); The Cambridge Dictionary of Science and Technology (Walker ed., 1988); The Glossary of Genetics, 5th Ed., R. Rieger et al. (eds.), Springer Verlag (1991); and Hale & Marham, The Harper Collins Dictionary of Biology (1991). As used herein, the following terms have the meanings ascribed to them below, unless specified otherwise.

[0012] As used herein, "adapter" refers to a nucleic acid sequence in an amplification primer that is complementary to the sequence of a nucleic acid molecule used to prime downstream sequencing reactions.

[0013] The term "allelic dropout" refers to the loss of one allele during amplification, resulting in apparent homozygosity. Nucleotide variation, cytosine methylation, and nucleic acid structure in the primer binding site of only one allele can cause allelic dropout when primer binding to the primer binding site is inhibited or reduced. For example, G-quadruplexes (secondary structures formed from stacks of G-quartets) present in the primer binding sites of an allele can prevent efficient priming of the template nucleic acid and lead to allelic dropout.

[0014] By "alternative allele" is meant an allele other than a reference allele. An alternative allele will have genetic variation that is not present in the reference allele. In some embodiments, a reference allele is a wildtype allele. A reference allele may differ between different populations, races, or ethnicities. Genetic variation present in an alternative allele can be nucleotide variation (i.e., a transition or a transversion), an insertion, or a deletion. An alternative allele may have a silent variant or mutation, a missense variant or mutation, or a nonsense variant or mutation.

[0015] By "alternative allele fraction" is meant the frequency of an allele, other than a reference allele, in a population of cells in an individual. The alternative allele

fraction is often less than that of the reference allele fraction, especially when the reference allele is a wildtype allele.

[0016] By "amplicon" is meant the product of an amplification reaction.

[0017] By "amplification bias" is meant a tendency for a nucleic acid amplification reaction to yield a particular amplicon. Amplification bias is often associated with inefficient primer binding. For example, if a primer's nucleic acid sequence is less complementary to the sequence of a template nucleic acid, the primer will be less likely to bind to the template than a primer having a more complementary sequence. Variants present in the primer binding site of a template nucleic acid may result in conformational or structural changes to the nucleic acid molecule that inhibit primer binding. Other variants or modifications (e.g., methylated nucleic acid residues) present in the primer binding site or elsewhere in the nucleic acid molecule can also cause to amplification bias. Amplification bias may result in underrepresentation of an allele or allelic dropout.

[0018] By "analog" is meant a molecule that is not identical, but has analogous functional or structural features to a naturally occurring molecule. For example, a polynucleotide analog retains the biological activity of a corresponding naturally-occurring polynucleotide while having certain modifications that enhance the analog's function relative to a naturally occurring polynucleotide. Such modifications could increase the polynucleotide's affinity for DNA, half-life, and/or nuclease resistance, an analog may include an unnatural nucleotide or amino acid.

[0019] By "bin" is meant a collection of sequencing reads that are substantially identical. In some instances, a bin comprises sequences reads that have the same index sequence or UMI sequence.

[0020] The phrase "biological sample" as used herein refers to a sample taken from a biological source and includes, but is not limited to, blood, serum, plasma, sputum, lavage fluid, cerebrospinal fluid, urine, semen, sweat, tears, tissue biopsy, and saliva. As used herein, the terms "blood," "plasma," and "serum" expressly encompass fractions or processed portions thereof.

[0021] In this disclosure, "comprises," "comprising," "containing" and "having" and the like can have the meaning ascribed to them in U.S. Patent law and can mean "includes," "including," and the like; "consisting essentially of" or "consists essentially" likewise has the meaning ascribed in U.S. Patent law and the term is open-ended, allowing for the presence of more than that which is recited so long as basic or novel characteristics of that which is recited is not changed by the presence of more than that which is recited, but excludes prior art embodiments.

[0022] By "demultiplex" is meant a process in which sequence reads generated from different amplicons are segregated into groups based on at least one characteristic unique to each group. For example, the index sequence of a primer can be used to segregate the sequence reads.


[0023] The term "denaturing," as contemplated herein, refers to removing impediments to primer binding from a nucleic acid. For example, denaturing includes removing conformational or structural properties of a nucleic acid or separating a nucleic acid duplex into single strands. Denaturing is facilitated by exposing the duplex to at least one denaturing condition or agent. Denaturing conditions are well known in the art. In one embodiment, a nucleic acid duplex is denatured by exposing it to a temperature that is

above the melting temperature (T_m) of the duplex. In certain embodiments, a nucleic acid may be denatured by exposing it to a temperature of at least 90° C. for a sufficient amount of time to denature the nucleic acid molecule. In some embodiments, a denaturing agent may include a chemical additive that facilitates denaturation, for example, sodium hydroxide or urea.

[0024] “Detect” refers to discovering or identifying the presence, absence, or amount of an analyte (e.g., genetic variation) to be detected.

[0025] By “detectable label” is meant a composition that when linked to a molecule of interest renders the latter detectable, via spectroscopic, photochemical, biochemical, immunochemical, or chemical means. For example, useful labels include radioactive isotopes, magnetic beads, metallic beads, colloidal particles, fluorescent dyes, electron-dense reagents, enzymes (for example, as commonly used in an ELISA), biotin, digoxigenin, or haptens.

[0026] “DMSO” refers to dimethyl sulfoxide, which has the following structure:

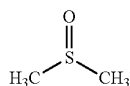

[0027] The term “enrich,” as used herein, refers to the process of further amplifying nucleic acid amplicons. In some embodiments, enrichment of nucleic acid amplicon allows for more efficient detection and quantifying of genetic variants having very low alternative allele frequency relative to detecting and quantifying genetic variants with very low alternative allele frequency in non-enriched nucleic acid amplicons.

[0028] By “GC buffer” is meant a reagent designed to optimize the ionic environment of an amplification reaction of a nucleic acid molecule having an enriched guanine/cytosine sequence.

[0029] “Germline allele” means an allele specific to germ cells or progenitors thereof.

[0030] “Hybridization” means hydrogen bonding, which may be Watson-Crick, Hoogsteen or reversed Hoogsteen hydrogen bonding, between complementary nucleobases. For example, adenine and thymine are complementary nucleobases that pair through the formation of hydrogen bonds.

[0031] By “index sequence” or “barcode” is meant a portion of a nucleic acid molecule that allows grouping or demultiplexing of sequencing reads. For example, an index sequence enables the segregation of sequence reads into bins, wherein each bin comprises sequence reads of amplicons generated from the primer pair having the index sequence. In some embodiments, each primer pair used in the presently disclosed methods has a unique index sequence.

[0032] As used herein, “interrogate” refers to obtaining nucleotide sequence information for a nucleic acid molecule.

[0033] The terms “isolated,” “purified,” or “biologically pure” refer to material that is free to varying degrees from components which normally accompany it as found in its native state. “Isolate” denotes a degree of separation from original source or surroundings. “Purify” denotes a degree

of separation that is higher than isolation. A “purified” or “biologically pure” nucleic acid is sufficiently free of other materials such that any impurities do not materially affect the biological properties of the nucleic acid or cause other adverse consequences. That is, a nucleic acid of this disclosure is purified if it is substantially free of cellular material, viral material, or culture medium. Purity and homogeneity are typically determined using analytical chemistry techniques, for example, polyacrylamide gel electrophoresis or high-performance liquid chromatography. The term “purified” can denote that a nucleic acid gives rise to essentially one band in an electrophoretic gel.

[0034] By “isolated polynucleotide” is meant a nucleic acid (e.g., a DNA) that is free of the genes which, in the naturally-occurring genome of the organism from which the nucleic acid molecule of the disclosure is derived, flank the gene. The term therefore includes, for example, a recombinant DNA that is incorporated into a vector; into an autonomously replicating plasmid or virus; or into the genomic DNA of a prokaryote or eukaryote; or that exists as a separate molecule (for example, a cDNA or a genomic or cDNA fragment produced by PCR or restriction endonuclease digestion) independent of other sequences. In addition, the term includes an RNA molecule that is transcribed from a DNA molecule, as well as a recombinant DNA that is part of a hybrid gene encoding additional polypeptide sequence.

[0035] “Isothermal” refers to a process incubated at about a constant temperature. For example, some isothermal amplification reactions are carried out at about 65° C. An isothermal temperature may depart from an intended temperature by not more than about 10% or 5° C., whichever is greater. An isothermal reaction may include an initial incubation at a higher temperature (“a hot start”). A hot start may comprise incubating the amplification reaction at a temperature sufficient to denature a region of interest on a nucleic acid molecule or to active a reagent (i.e., a polymerase).

[0036] By “marker” is meant any protein or polynucleotide associated with a disease or disorder.

[0037] As used herein, “mosaic” refers to two or more cells or populations of cells with different genotypes within an individual subject. For example, “somatic mosaicism” refers to two or more genotypically distinct somatic cells or populations of somatic cells in an individual. “Germline mosaicism” occurs when two or more genotypically distinct germ cells or populations of germ cells are present in an individual. Germline mosaicism generally arises after a mutation gives rise to a genotypically distinct gamete.

[0038] The term “Next Generation Sequencing (NGS)” refers to massive parallel sequencing of clonally amplified molecules or single nucleic acid molecules. “Massive parallel sequencing” refers to simultaneously performing more than 1000 separate, parallel sequencing reactions. Non-limiting examples of NGS include sequencing-by-synthesis using reversible dye terminators, sequencing-by-ligation, and electronic detection sequencing methods. Electronic detection sequencing methods include those used in the Ion Torrent sequencing strategy (ThermoFisher Scientific) or MiSeq platform (Illumina), wherein changes in pH are detected when a nucleotide is incorporated into a nucleic acid strand resulting in release of a hydrogen ion.

[0039] The terms “nucleic acid” and “nucleic acid molecule,” are used interchangeably herein and refer to a compound comprising a nucleobase and an acidic moiety,

e.g., a nucleoside, a nucleotide, or a polymer of nucleotides. Typically, polymeric nucleic acids, e.g., nucleic acid molecules comprising three or more nucleotides are linear molecules, in which adjacent nucleotides are linked to each other via a phosphodiester linkage. In some embodiments, “nucleic acid” refers to individual nucleic acid residues (e.g. nucleotides and/or nucleosides). In some embodiments, “nucleic acid” refers to an oligonucleotide chain comprising three or more individual nucleotide residues. As used herein, the terms “oligonucleotide” and “polynucleotide” can be used interchangeably to refer to a polymer of nucleotides (e.g., a string of at least three nucleotides). In some embodiments, “nucleic acid” encompasses RNA as well as single and/or double-stranded DNA. Nucleic acids may be naturally occurring, for example, in the context of a genome, a transcript, an mRNA, tRNA, rRNA, siRNA, snRNA, a plasmid, cosmid, chromosome, chromatid, or other naturally occurring nucleic acid molecule. On the other hand, a nucleic acid molecule may be a non-naturally occurring molecule, e.g., a recombinant DNA or RNA, an artificial chromosome, an engineered genome, or fragment thereof, or a synthetic DNA, RNA, DNA/RNA hybrid, or including non-naturally occurring nucleotides or nucleosides. Furthermore, the terms “nucleic acid,” “DNA,” “RNA,” and/or similar terms include nucleic acid analogs, e.g., analogs having other than a phosphodiester backbone. Nucleic acids can be purified from natural sources, produced using recombinant expression systems and optionally purified, chemically synthesized, etc. Where appropriate, e.g., in the case of chemically synthesized molecules, nucleic acids can comprise nucleoside analogs such as analogs having chemically modified bases or sugars, and backbone modifications. A nucleic acid sequence is presented in the 5' to 3' direction unless otherwise indicated. In some embodiments, a nucleic acid is or comprises natural nucleosides (e.g. adenosine, thymidine, guanosine, cytidine, uridine, deoxyadenosine, deoxythymidine, deoxyguanosine, and deoxycytidine); nucleoside analogs (e.g., 2-aminoadenosine, 2-thiothymidine, inosine, pyrrolo-pyrimidine, 3-methyl adenosine, 5-methylcytidine, 2-aminoadenosine, C5-bromouridine, C5-fluorouridine, C5-iodouridine, C5-propynyl-uridine, C5-propynyl-cytidine, C5-methylcytidine, 2-aminoadenosine, 7-deazaadenosine, 7-deazaguanosine, 8-oxoadenosine, 8-oxoguanosine, 0(6)-methylguanine, and 2-thiocytidine); chemically modified bases; biologically modified bases (e.g., methylated bases); intercalated bases; modified sugars (e.g., 2'-fluororibose, ribose, 2'-deoxyribose, arabinose, and hexose); and/or modified phosphate groups (e.g., phosphorothioates and 5'-N-phosphoramidite linkages).

[0040] Nucleic acid molecules assayed using the methods described herein need not be 100% identical with an endogenous nucleic acid sequence, but will typically exhibit substantial identity. Polynucleotides having “substantial identity” to an endogenous sequence are typically capable of hybridizing with at least one strand of a double-stranded nucleic acid molecule. Nucleic acid molecules useful in the methods of the disclosure include any nucleic acid molecule that encodes a polypeptide of the disclosure or a fragment thereof. Such nucleic acid molecules need not be 100% identical with an endogenous nucleic acid sequence, but will typically exhibit substantial identity. Polynucleotides having “substantial identity” to an endogenous sequence are typically capable of hybridizing with at least one strand of a double-stranded nucleic acid molecule. By “hybridize” is

meant pair to form a double-stranded molecule between complementary polynucleotide sequences (e.g., a gene described herein), or portions thereof, under various conditions of stringency. (See, e.g., Wahl, G. M. and S. L. Berger (1987) *Methods Enzymol.* 152:399; Kimmel, A. R. (1987) *Methods Enzymol.* 152:507).

[0041] For example, stringent salt concentration will ordinarily be less than about 750 mM NaCl and 75 mM trisodium citrate, less than about 500 mM NaCl and 50 mM trisodium citrate, or about 250 mM NaCl and 25 mM trisodium citrate. Low stringency hybridization can be obtained in the absence of organic solvent, e.g., formamide, while high stringency hybridization can be obtained in the presence of at least about 35% formamide, and in some embodiments, at least about 50% formamide. Stringent temperature conditions will ordinarily include temperatures of at least about 30° C. at least about 37° C., or at least about 42° C. Varying additional parameters, such as hybridization time, the concentration of detergent, e.g., sodium dodecyl sulfate (SDS), and the inclusion or exclusion of carrier DNA, are well known to those skilled in the art. Various levels of stringency are accomplished by combining these various conditions as needed. In one embodiment, hybridization will occur at 30° C. in 750 mM NaCl, 75 mM trisodium citrate, and 1% SDS. In another embodiment, hybridization will occur at 37° C. in 500 mM NaCl, 50 mM trisodium citrate, 1% SDS, 35% formamide, and 100 µg/ml denatured salmon sperm DNA (ssDNA). In yet another embodiment, hybridization will occur at 42° C. in 250 mM NaCl, 25 mM trisodium citrate, 1% SDS, 50% formamide, and 200 µg/ml ssDNA. Useful variations on these conditions will be readily apparent to those skilled in the art.

[0042] For most applications, washing steps that follow hybridization will also vary in stringency. Wash stringency conditions can be defined by salt concentration and by temperature. As above, wash stringency can be increased by decreasing salt concentration or by increasing temperature. For example, stringent salt concentration for the wash steps will comprise less than about 30 mM NaCl and 3 mM trisodium citrate or less than about 15 mM NaCl and 1.5 mM trisodium citrate. Stringent temperature conditions for the wash steps will ordinarily include a temperature of at least about 25° C., at least about 42° C., or at least about 68° C. In some embodiments, wash steps will occur at 25° C. in 30 mM NaCl, 3 mM trisodium citrate, and 0.1% SDS. In other embodiments, wash steps will occur at 42° C. in 15 mM NaCl, 1.5 mM trisodium citrate, and 0.1% SDS. In other embodiments, wash steps will occur at 68° C. in 15 mM NaCl, 1.5 mM trisodium citrate, and 0.1% SDS. Additional variations on these conditions will be readily apparent to those skilled in the art. Hybridization techniques are well known to those skilled in the art and are described, for example, in Benton and Davis (*Science* 196:180, 1977); Grunstein and Hogness (*Proc. Natl. Acad. Sci., USA* 72:3961, 1975); Ausubel et al. (*Current Protocols in Molecular Biology*, Wiley Interscience, New York, 2001); Berger and Kimmel (*Guide to Molecular Cloning Techniques*, 1987, Academic Press, New York); and Sambrook et al., *Molecular Cloning: A Laboratory Manual*, Cold Spring Harbor Laboratory Press, New York.

[0043] As used herein, “obtaining” as in “obtaining an agent” includes synthesizing, purchasing, or otherwise acquiring the agent.

[0044] By “overlapping amplicons” is meant two or more amplicons that comprise a shared nucleic acid sequence but have at least one different terminal sequence.

[0045] “Polymerase” refers to an enzyme capable of catalyzing nucleic acid synthesis. A polymerase can be a DNA polymerase or an RNA polymerase. A polymerase can be characterized by its error rate, or the rate at which the polymerase inserts an incorrect nucleotide into the nucleic acid molecule it is synthesizing. In some embodiments, a polymerase can be a high-fidelity polymerase, which has a much lower error rate than a reference polymerase. A non-limiting example of a reference polymerase is Taq polymerase.

[0046] “Pooling,” as used herein, means combining multiple amplification reactions or groups of reactions. Pooling is synonymous with multiplexing.

[0047] By “portion” is meant a segment of an intact nucleic acid molecule. This portion contains, in some embodiments, at least 10%, 20%, 30%, 40%, 50%, 60%, 70%, 80%, or 90% of the entire length of the reference nucleic acid molecule. A portion may contain 10, 20, 30, 40, 50, 60, 70, 80, 90, or 100, 200, 300, 400, 500, 600, 700, 800, 900, or 1000 nucleotides.

[0048] The term “read,” “sequence read,” or “sequencing read” refers to sequencing data from a region of a nucleic acid molecule obtained from a single nucleic acid molecule. A read represents a short sequence of contiguous bases in the nucleic acid molecule and may be depicted, for example, as a chromatogram or as a linear string of letters that represent the nitrogenous bases of the nucleotide sequence, wherein A=adenine; G=guanine; C=cytosine; T=thymine; U=uracil; R=purine (A or G); Y=pyrimidine (C or T); N=any nucleotide; W=A or T; S=G or C; K=G or T; B=Not A; H=Not G; D=Not C; and V=Not T.

[0049] “Reduces” or “increases” refers to a negative or positive alteration, respectively, of at least 10%, 25%, 50%, 75%, or 100%.

[0050] By “reference” is meant a standard or control condition.

[0051] A “reference sequence” is a defined sequence used for sequence comparison. A reference sequence may be a subset of or the entirety of a specified sequence; for example, a segment of a full-length gene sequence, or the complete gene sequence. For nucleic acids, the length of the reference nucleic acid sequence will generally be at least about 50 nucleotides, at least about 60 nucleotides, at least about 75 nucleotides, about 100 nucleotides, or even about 300, 400, or 500 nucleotides or any integer thereabout or therebetween. In some embodiments, the length of the reference nucleic acid sequence will be less than 50 nucleotides. In some embodiments, the reference nucleic acid sequence will be more than 500 nucleotides.

[0052] The term “sequence variant,” as used herein, refers to an alteration in a sequence relative to a reference sequence. In one embodiment, a nucleotide sequence variant comprises one or more alterations relative to a reference nucleotide sequence. In some embodiments, the reference sequence is a consensus sequence. Optimally aligned sequencing reads obtained from multiple individuals of the same species or a population thereof, or multiple sequencing reads for the same individual, may be used to produce a consensus sequence. As contemplated herein, a “consensus sequence” refers to a nucleotide sequence that comprises the

base most in common among all the sequencing reads at each nucleotide in the sequence.

[0053] In some embodiments, a sequence variant represents a variation relative to corresponding sequences in the same sample. In some embodiments, the sequence variant occurs with a low frequency (i.e., at least <1%) in the population (also referred to as a “rare variant”). For example, the sequence variant may occur with a frequency of about or less than about 5%, 4%, 3%, 2%, 1.5%, 1%, 0.75%, 0.5%, 0.25%, 0.1%, 0.075%, 0.05%, 0.04%, 0.03%, 0.02%, 0.01%, 0.005%, 0.001%, or lower. In some embodiments, the sequence variant occurs with a frequency above about 0.1%. In some embodiments, the sequence variant occurs at a frequency of above about 0.0025%.

[0054] By “somatic allele” is meant an allele specific to a non-germline cell (i.e., somatic cell).

[0055] By “somatic event” is meant the acquisition of a genetic variant by a somatic cell.

[0056] By “subject” is meant a mammal, including a human or a non-human mammal, such as a bovine, equine, canine, ovine, feline, or rodent (e.g., mouse, rat).

[0057] By “substantially identical” is meant a polypeptide or nucleic acid molecule exhibiting at least 50% identity to a reference amino acid sequence (for example, any one of the amino acid sequences described herein) or nucleic acid sequence (for example, any one of the nucleic acid sequences described herein). In some embodiments, such a sequence is at least 60%, 80% or 85%, 90%, 95% or even 99% identical at the amino acid level or nucleic acid to the sequence used for comparison.

[0058] Sequence identity is typically measured using sequence analysis software (for example, Sequence Analysis Software Package of the Genetics Computer Group, University of Wisconsin Biotechnology Center, 1710 University Avenue, Madison, Wis. 53705, BLAST, BESTFIT, GAP, or PILEUP/PRETTYBOX programs). Such software matches identical or similar sequences by assigning degrees of homology to various substitutions, deletions, and/or other modifications. In an exemplary approach to determining the degree of identity, a BLAST program may be used, with a probability score between e^{-3} and e^{-100} indicating a closely related sequence.

[0059] The term “tissue” refers to a group or layer of similarly specialized cells, which together perform certain special functions. The term “tissue-specific” refers to a source or defining characteristic of cells from a specific tissue.

[0060] By “unique molecular identifier (UMI)” is meant a distinct nucleic acid sequence that individualizes each primer used in an amplification reaction. For example, 500 primers having identical complementary nucleic acid sequences will have 500 different UMIs. UMIs facilitate the detection and removal of redundant sequencing reads.

[0061] Unless specifically stated or obvious from context, as used herein, the term “or” is understood to be inclusive. Unless specifically stated or obvious from context, as used herein, the terms “a,” “an,” and “the” are understood to be singular or plural.

[0062] Unless specifically stated or obvious from context, as used herein, the term “about” is understood as within a range of normal tolerance in the art, for example within 2 standard deviations of the mean. About can be understood as within 10%, 9%, 8%, 7%, 6%, 5%, 4%, 3%, 2%, 1%, 0.5%,

0.1%, 0.05%, or 0.01% of the stated value. Unless otherwise clear from context, all numerical values provided herein are modified by the term about.

[0063] Ranges provided herein are understood to be shorthand for all the values within the range. For example, a range of 1 to 50 is understood to include any number, combination of numbers, or sub-range from the group consisting 1, 2, 3, 4, 5, 6, 7, 8, 9, 10, 11, 12, 13, 14, 15, 16, 17, 18, 19, 20, 21, 22, 23, 24, 25, 26, 27, 28, 29, 30, 31, 32, 33, 34, 35, 36, 37, 38, 39, 40, 41, 42, 43, 44, 45, 46, 47, 48, 49, or 50.

[0064] The recitation of a listing of chemical groups in any definition of a variable herein includes definitions of that variable as any single group or combination of listed groups. The recitation of an embodiment for a variable or aspect herein includes that embodiment as any single embodiment or in combination with any other embodiments or portions thereof.

[0065] Any compositions or methods provided herein can be combined with one or more of any of the other compositions and methods provided herein.

BRIEF DESCRIPTION OF THE DRAWINGS

[0066] FIGS. 1A to 1C are schematic diagrams illustrating the primer design strategy used in the presently disclosed methods. FIG. 1A is a schematic diagram illustrating overlapping amplicons that provide redundant coverage of a variant of interest (A/G). Primer 1, Primer 2, and Primer 3 refer to the pairs of forward and reverse primers (depicted at the termini of the intervening line). The intervening line represents the nucleic acid sequence to be amplified. “SNV” refers to single nucleotide variant. FIG. 1B is a schematic diagram of three amplicons, wherein “Adapter 1” and “Adapter 2” refer to the adapter sequences upstream from the primer’s complementary nucleotide sequence (“Forward” or “Reverse”). Each reverse primer has one of three index sequences. FIG. 1C is a schematic diagram of three amplicons that comprise a unique molecular identifier (UMI).

[0067] FIG. 2 comprises three panels of aligned sequencing reads, wherein each panel comprises sequencing reads of amplicons generated from one of three amplification reactions. The top and bottom panels each show alternate allele fractions of a detected variant of approximately 50%. The middle panel shows an alternate allele fraction of only 3%, which indicates allelic dropout.

[0068] FIG. 3 is an illustration of capturing and enriching amplified nucleic acids.

[0069] FIG. 4 is a schematic diagram of a method for detecting low frequency variants in a nucleic acid molecule. Throughout the figures, QC denotes quality control and AAF denotes alternative allele fractions.

[0070] FIG. 5A is a schematic diagram of a method for detecting and characterizing low frequency variants. CI denotes confidence interval. FIG. 5B is a diagram illustrating an optional quality control step that can be added to the method depicted in FIG. 5A.

[0071] FIG. 6 is a chart summarizing an Ion Torrent Next Generation Sequencing run and the data generated therefrom.

[0072] FIG. 7 is an illustration of demultiplexing sequencing data.

[0073] FIG. 8 is data output illustrating sequencing errors generated using the Ion Torrent platform. Specifically, the

data presented illustrates how sequencing errors (i.e., indels) are processed using the disclosed methods.

[0074] FIG. 9 is an illustration of sequencing reads, wherein the ends of each read (i.e., the primer sequences) are easily observed.

[0075] FIG. 10A is an illustration of the reproducibility observed in aligned sequencing data of a germline event. The illustration depicts three panels of aligned sequence data indicating the presence of a variant at base pair number 14,234,400. FIG. 10B is an illustration of a detected mutation.

[0076] FIGS. 11A to 11G graphically illustrate quality control assessment of amplification products generated using the methods as described herein. FIG. 11A is an electronically generated gel image of products of an amplification reaction performed according to the methods described herein. Lane (L) 1 comprises a control sample “Control-6-U” that was not amplified using the methods disclosed herein. Lane 2 comprises amplification products generated using a single amplification (20 cycles) protocol as described herein. Lane 3 comprises amplification products using a two-amplification protocol (first amplification=8 cycles; second amplification=20 cycles). “Bio” indicates the first-round amplification products were biotinylated. Lane 4 comprises amplification products generated using a two-amplification protocol (first amplification=10 cycles; second amplification=20 cycles). Lane 5 comprises amplification products generated using a two-amplification protocol (first amplification=10 cycles; second amplification=20 cycles). “Amp” indicates the first-round reaction products were not biotinylated. “[s]” refers to seconds. FIG. 11B is a graph illustrating the fluorescent peaks detected when analyzing the control reaction “Control-6-U” using the Bioanalyser 2100. FIG. 11C is a graph illustrating the fluorescent peaks detected when analyzing the “20X-Norm” reaction using the Bioanalyser 2100. FIG. 11D is a graph illustrating the fluorescent peaks detected when analyzing the “8X_20X_Bio” reaction using the Bioanalyser 2100. FIG. 11E is a graph illustrating the fluorescent peaks detected when analyzing the “10X_20X_Bio” reaction using the Bioanalyser 2100. FIG. 11F is a graph illustrating the fluorescent peaks detected when analyzing the “10X_20X_Amp” reaction using the Bioanalyser 2100. FIG. 11G is a graph illustrating the fluorescent peaks detected when analyzing the “Exo 8X_20X_RD” reaction using the Bioanalyser 2100. This reaction was purified using the ExoSAP protocol described herein after amplifying a target nucleic acid using a two-amplification protocol as used herein. In this sample, the target nucleic acid was amplified with a first reaction comprising 8 cycles and then a subsequent amplification reaction comprising 20 cycles.

[0077] FIG. 12 is a graph depicting a TapeStation analyzer’s quality control assessment of the products generated in an amplification reaction. The “upper” and “lower” peaks are the control peaks, and the “283” peak represents the amplification reaction products.

[0078] FIG. 13 is a graph illustrating the accuracy and reproducibility of the present methods to detect variants and provide accurate alternative allele fractions.

[0079] FIG. 14 is a graph illustrating the accuracy and reproducibility of the present methods to detect low frequency variants and provide accurate alternative allele fractions (i.e., AAF<1%).

[0080] FIG. 15 is a graph of a deleterious missense mosaic variant detected in the CACNA1A gene of a single individual.

[0081] FIG. 16 is a graph of number of germline heterozygous single nucleotide having a particular variant (alternate allele fraction (VAF)).

[0082] FIGS. 17A to 17D are graphs and figures explaining asymmetric cell contribution. FIG. 17A is a graph showing asymmetrical cell contributions to brain development during early embryonic development. FIG. 17B is an illustration of the different branches of early phylogeny at which mutations may be acquired. FIG. 17C is a graph showing poor stability of the asymmetric parameter α_1 estimated from the 2nd cell generation compared to only one asymmetric cell division. FIG. 17D is a graph showing the confidence interval for the asymmetric cell contribution parameter.

[0083] FIGS. 18A-18D illustrate that the presently described methods accurately measure AAFs as low as 0.01% when using a 50 ng of genomic DNA. FIG. 18A is a graph showing the correlation of expected and measured AAFs up to 60% for samples comprising 50 ng of DNA. FIG. 18B is a graph showing the correlation of expected and measured AAFs between 0 and 1.0%. FIG. 18C is a graph showing the correlation of expected and measured AAFs up to 60% for samples comprising 25 ng of DNA. FIG. 18D is a graph showing the correlation of expected and measured AAFs between 0 and 1.0%.

[0084] FIG. 19A is a graph correlating the AAF's of single nucleotide variants determined using whole genome sequencing (WGS) and triple-primer PCR sequencing (Trip-Seq). FIG. 19B is a graph correlating the AAF's of indels determined using whole genome sequencing (WGS) and triple-primer PCR sequencing (Trip-Seq). FIG. 19C is a graph showing the correlation of expected and measured AAFs when consistent AAFs are required across multiple unique primer sets. FIG. 19D is a graph of the expected and measured AAFs when triple-primer PCR sequencing is applied to a large set of tissues derived DNA samples for detections of novel mutations in a given gene.

DETAILED DESCRIPTION OF THE INVENTION

[0085] The present disclosure features methods for detecting and quantifying genetic variants in a sample.

[0086] The invention is based, at least in part, on the discovery of triple primer PCR sequencing ("TriPP-seq"), which provides a highly sensitive, low-cost approach for detecting and validating mutation on a highly scalable system. Mosaic mutations in somatic or germline cells contribute to a wide range of human disorders. As such, their identification and accurate allelic fraction quantification from tissue-derived and cell-free DNA are essential for clinical diagnoses and early detection of cancers. However, rapid, low-cost detection and validation of ultra-low alternate allelic fraction (AAF) mutations has traditionally required expensive and low throughput methods that have limited widespread testing. Recent methods, (e.g., ddPCR) have shown great promise for detection and validating known mutations at very low AAFs, but remain low-throughput due to allele-specific optimization.

[0087] Accordingly, the present disclosure features methods for detecting low frequency genetic variation. The present disclosure's novel approach is based on generating

deep coverage of overlapping amplicons of a target nucleic acid sequence. Because the primers used in the reactions are designed to allow discernment and segregation of the overlapping amplicons, the sequencing data can be segregated into groups, and analysis of the sequencing data can be performed in parallel. The methods provide not only deep coverage of the target nucleic acid, but also a cost-effective means of characterizing and validating sequencing results.

[0088] Recently, the important roles of somatic mutations beyond cancer are becoming more appreciated with discoveries of somatic mutations across a wide range of neurodevelopmental, overgrowth, and hematological disorders. Even more, the presence of somatic mutations in healthy cells and individuals are associated with normal development and aging and are, therefore, a powerful tool for understanding how cells divide and form complex organs like the human brain. Finally, with the detection of cell-free DNA (e.g., fetal and tumor), it is becoming possible for early detection of disease, tracking of disease recurrence in cancers, and even non-invasive prenatal genetic testing where mutations of the placenta are detected in the pregnant mother's blood sample. The rapid advancements in sequencing technologies and interest in genetic mutation present at low alternate allelic fraction (i.e., ratio of DNA fragments carrying the mutation to those with the wild-type allele in a given sample; AAF) poses some major challenges for both the clinical and research communities related to the sensitivity to detect mutations, false positives, and the precision of the assessed AAFs. These challenges are often confounded by the inability to directly assess tissues with the highest AAFs, as is the case with brain tissue, or by limited or degraded DNA samples, as is typical for cell free DNA.

[0089] While germline mutations are relatively easy to detect with small amounts of DNA with variable qualities using WES, WGS, targeted gene panels, and traditional Sanger sequencing due to the equal fractions of mutant to wild-type alleles (50% AAF) in a given DNA sample, the AAF of a somatic mutation will depend on the given tissue, cell type, and the stage in development at which the mutation arose. Traditional WGS and WES sequencing in both the research and clinical diagnostic settings are optimized to identify germline events, but often lack the sequencing depth to robustly detect low-AAF variants. However, many recently improvements allow for robust detection of mutations present at greater than 0.1% AAF. These tools often employ strategies such as molecular barcoding, increased read depth, and reduced use of PCR to mitigate sequencing-induced errors while improving sensitivity. Despite these measures, the identification of somatic alleles, particularly those at very low AAFs has an elevated false positive rate compared to germline mutations. Therefore, while essential, the validation of large numbers of somatic alleles is often challenging due to many factors like assay costs, throughput, and sensitivity limitations.

[0090] The methodology utilized to accurately detect or validate somatic mutations have rapidly advanced in the last few years. The challenge for validating or measuring low AAFs is multifaceted, spanning sequencing platforms, inherent error rates of polymerases, and locus specific challenges. Each of these result in additional errors and skewing of AAFs, which can mask or alter the detected AAF in each assay. The utilization of PCR to amplify the genomic loci without inducing additional mutations and maintain the original AAFs has been improved using improved poly-

merases with proofreading capabilities and, in some cases, unique molecular barcodes for each DNA fragment. Additionally, errors can occur during sequencing on both the Illumina and Ion Torrent platforms. For example, in one study, the Ion Torrent had an error rate ~0.05% for SNVs but ~1.5% for indels while the on the Illumina MiSeq had 0.1% errors for SNVs and 0.7% for indels.

[0091] The original methods used employed either pyrosequencing or bacterial cloning followed by sanger sequencing of hundreds or thousands of individual bacterial colonies to measure a single mutation. These methods, while accurate and robust, were often cost-prohibitive, less scalable to large numbers of mutations, and were less sensitive for mutations below 5% AAF. These methods were recently succeeded by the advancement of digital droplet PCR, ddPCR, where an allele-specific PCR conditions are designed to allow for the measurement of mutation positive and negative DNA fragments in thousands of droplets. This method is routinely considered a gold standard for validation of somatic alleles in both research and clinical settings, but each allele requires the development of a custom assay, validation and optimization prior to use. The ddPCR assay can accurately detect AAFs below 0.5%, but its sensitivity relies on the quantity and concentration of input DNA and the number of positive droplets formed in each reaction. Despite its great success, the use of ddPCR is somewhat limited as it remains limited by scalability, the potential for allelic dropout, and the ability to design allele-specific primers, which is more challenging in repetitive regions and for small indels.

[0092] The growing consensus that somatic mutations might underlie a wide range of clinical phenotypes ranging from cancer risk to severe neurodevelopmental and overgrowth conditions suggests that a robust method for both detection and validation of alleles and their mosaic fraction in the body is essential. Here, an improved strategy that aims to mitigate the previously stated limitations for assessing somatic mutations is presented. This strategy, which can be referred to as triple-primer PCR, relies on the power of designing and running at least 3 unique, nonoverlapping amplicons over a suspected mutation. Through independently analyzing each amplicon, the impact of allelic dropout, amplification bias, sequencing and PCR induced artifacts, and general optimization challenges, are markedly reduced while achieving the highest sensitivity to accurately detect ultra-low allelic fractions below 0.1% regardless of tissue origin. As described, below, this triple-primer PCR sequencing method allows for additional improvements to future improve accuracy through incorporations of molecular barcoding and improved purification processes.

Primers

[0093] Nucleic acid amplification according to the presently disclosed methods requires at least two pairs of primers and in some embodiments, at least three pairs of primers. Each pair of primers comprises a forward and a reverse primer, and each primer comprises a complementary nucleic acid sequence that is at least 85% complementary to a nucleic acid sequence (i.e., the primer binding site) on a template nucleic acid molecule. The primers of each pair define the termini of an amplicon that is generated by an amplification reaction, and the region of the amplicon between the termini comprises the target nucleic acid sequence. The combined length of the primers and the target sequence is referred to as the amplicon length. Amplicon

length is typically between about 150 and about 500 nucleotides. In some embodiments, the length of the amplicon is about 150, 200, 250, 300, 350, 400, 450, 500, or any integer in-between, nucleotides. In some embodiments, the length of the amplicon is less than 150 nucleotides. In some embodiments, the length of the amplicon is greater than 500 nucleotides. Each primer has a unique nucleic acid sequence that can bind to a complementary primer binding site on the template nucleic acid.

[0094] Amplicons generated by amplification reactions using one of the primer pairs will be distinguishable from other amplicons generated by amplification reactions that use different primer pairs due to the length and sequence of the amplicon (FIG. 1A). Each amplicon will include the target nucleic acid sequence, and because the primers are designed to generate overlapping amplicons, each amplicon is at least partially redundant to the other amplicons. In other embodiments, only one primer of each pair will have a unique complementary nucleic acid sequence, such that the amplicons have either the same 5' terminus nucleic acid sequence and differing 3' terminus nucleic acid sequences or differing 5' terminus nucleic acid sequences and the same 3' terminus nucleic acid sequence.

[0095] A primer binding site in a template nucleic acid sequence may harbor a variant that impairs primer binding, which results in decreased amplification of the template harboring the variant and a loss of sequencing coverage of the allele. The resulting loss of coverage of a particular variant is allelic dropout. Referring to FIG. 2, three panels of sequencing data (derived from three sets of overlapping amplicons) show allelic dropout in the middle panel. To minimize allelic dropout in amplification reactions comprising one of three (or more) pairs of primers, at least two forward primers and at least two of the three reverse primers have different complementary nucleic acid sequences. If only two pairs of primers are used, both forward primers and both reverse primers should have unique complementary nucleic acid sequences.

[0096] In some embodiments, the complementary nucleic acid sequence of a primer is about 15, 16, 17, 18, 19, 20, 25, 30, 35, or even 40 nucleotides long. In some embodiments, the complementary nucleic acid sequence of a primer is between about 85% and about 100% complementary to a nucleic acid sequence in the template nucleic acid molecule. In some embodiments, the complementary nucleic acid sequence of the primer is between about 85%, 86%, 87%, 88%, 89%, 90%, 91%, 92%, 93%, 94%, 95%, 96%, 97%, 98%, 99%, and 100% complementary to a nucleic acid sequence in the template nucleic acid molecule. In some embodiments, wherein the complementary nucleic acid sequence of the primer is less than 100% complementary with a primer binding site in the template nucleic acid molecule, the mismatch nucleotide or nucleotides in the primer reside at least three bases from the 3' terminus of the primer. This allows for efficient binding at the terminus of the primer to the template molecule, which facilitates polymerase binding to the primer:template hybrid and extending the primer.

[0097] In some embodiments, a primer is comprised of DNA or RNA nucleotides. In some embodiments, a primer comprises at least one modified base. A modified base includes, but is not limited to, those nucleotide analogs described herein or a labeled nucleotide. In some embodiments, a primer may have a modified backbone comprising

at least one phosphorothioate linkage. In some embodiments, the primer comprises a label, such as, but not limited to, a fluorescent label, a radiolabel, a nanoparticle label, and/or a biotin label.

[0098] In some embodiments, each primer will have an adapter upstream from the complementary nucleic acid sequence. The adapter has a nucleic acid sequence that is complementary to a sequence of a nucleic acid molecule used in a downstream sequencing reaction. For example, the adapters used in some embodiments are designed to be compatible with Next Generation Sequencing including, but not limited to, Ion Torrent and MiSeq platforms. In some embodiments, the length of the adapter is between 8 and 20 nucleotides. In some embodiments, the length of the adapter is 8, 9, 10, 11, 12, 13, 14, 15, 16, 17, 18, 19, or 20 nucleotides. The adapter's sequence is designed to reduce or eliminate nonspecific binding of the adapter to the template nucleic acid molecule. In some embodiments, the adapter is designed to have a sequence that is not substantially complementary to any nucleic acid sequence present in the template nucleic acid molecule. In some embodiments, the adapter is designed to diverge from perfect complementarity with the template by 2, 3, or 4 or more nucleotides.

[0099] At least one primer in each pair also has an index sequence, or barcode (FIG. 1B). The index sequence allows for rapid identification of sequencing data generated from similar amplicons. The index sequence as contemplated herein can be between 8 and 30 nucleotides in length. For example, the index sequence contemplated herein may comprise 8, 9, 10, 11, 12, 13, 14, 15, 16, 17, 18, 19, 20, 21, 22, 23, 24, 25, 26, 27, 28, 29, or 30 nucleotides. Similar to the adapter, the index sequence is designed to reduce or eliminate nonspecific binding of it to the template nucleic acid molecule. In some embodiments, the index sequence comprises a nucleic acid sequence that is not substantially complementary to any nucleic acid sequence present in the template nucleic acid molecule. In some embodiments, the index sequence is designed to diverge from perfect complementarity with a nucleic acid sequence in the template nucleic acid molecule by 2, 3, or 4 or more nucleotides. In some embodiments, the index sequence is designed so that the most complementary sequence in the template has a conformation or structure that disfavors index sequence binding.

[0100] In some embodiments, at least one primer in each pair comprises a unique molecular identifier (UMI) (FIG. 1C). A UMI may allow for the detection of redundant sequencing reads. As contemplated herein, the UMI will comprise between 5 and 20 nucleotides. For example, the UMI contemplated herein may comprise 5, 6, 7, 8, 9, 10, 11, 12, 13, 14, 15, 16, 17, 18, 19, or 20 nucleotides. In some embodiments, no two primers will have the same UMI. Similar to the adapter and the index sequence, UMIs are designed to reduce or eliminate nonspecific binding of the UMIs to the template nucleic acid molecule. In some embodiments, the UMI comprises a nucleic acid sequence that is not substantially complementary to any nucleic acid sequence present in the template nucleic acid molecule. In some embodiments, the UMI is designed to diverge from perfect complementarity with the template by 2, 3, or 4 or more nucleotides. In some embodiments, the UMIs are designed so that the most complementary sequences in the template nucleic acid have a conformation that disfavors UMI binding.

[0101] There are approximately 1,000 possible sequences for a 5-nucleotide UMI, approximately 65,000 possible sequences for an 8-nucleotide UMI, approximately 1×10^6 possibilities for a 10-nucleotide UMI, and approximately 1×10^{12} possibilities for a 20-nucleotide UMI. Even if some UMIs are not suitable for the reasons given above, large UMI libraries can be produced for use in the presently disclosed methods. Use of nucleotide analogs increases the number of possible sequences for a UMI.

[0102] Table 1 characterizes five primer pairs used in the disclosed methods. In this table, "Chr. No." means chromosome number; "Ref" refers to the reference nucleotide; and "Alt" refers to the alternate nucleotide. Each of the primer pairs is designed to amplify a region containing a single nucleotide variant (the "allele start" and "allele end" are the same locus number). Three of the primer pairs on Table 1 (X:153579431-153579431/T/C-F1; X:153579431-153579431/T/C-F2; and X:153579431-153579431/T/C-F3) are used to interrogate a single nucleotide variant in the Filamin A (FLNA) gene on the X chromosome. The remaining two primer pairs (X:153579431-153579431/T/C-F1 and X:153579431-153579431/T/C-F2) are used to interrogate a single nucleotide variant in the SR-Related CTD Associated Factor 11 (SCAF-11) gene on chromosome 12. The amplicons generated in amplification reactions comprising the primer pairs disclosed in Table 1 will be about 220 to 260 nucleotides in length.

TABLE 1

PrimerID	Chr. No.	Allele Start	Allele End	Allele	Ref	Alt	Gene	Sample ID	Prod. Start	Prod. End	Insert Start	Insert End
1	X:153579431- 153579431/T/ C-F1	15357 9431	15357 9431	T	C		FLNA	PH4201	153579266	153579517	153579284	153579499
2	X:153579431- 153579431/T/ C-F2	15357 9431	15357 9431	T	C		FLNA	PH4201	153579289	153579555	153579311	153579536
3	X:153579431- 153579431/T/ C-F3	15357 9431	15357 9431	T	C		FLNA	PH4201	153579379	153579637	153579397	153579619
4	12:46321441- 46321441/T/ G-F1	46321 441	46321 441	T	G		SCAF11	PH4201	46321317	46321542	46321343	46321517
5	12:46321441- 46321441/T/ G-F2	46321 441	46321 441	T	G		SCAF11	PH4201	46321246	46321470	46321271	46321448

Primer ID	Forward	Reverse	Barcode No.	Primer barcode	type	Forward	UMI
1	X:153579431- 153579431/T/ C-F1	CAGGGCTCACC TTGGTC	ttaacggagcgc GCCA	1	Bar- code	CAAGGT GAGGCC CTG	No
2	X:153579431- 153579431/T/ C-F2	CTGTGACATAGC ACTCCTCCAG	tccggcttactGCAAAATC AGTGGCTCTCC	2	Bar- code	AGTGCT ATGTCAC AG	No
3	X:153579431- 153579431/T/ C-F3	AGGCTGGCTGGT TGCACCT	tctcattcagCTCCCTTCC TGCCACCTG	3	Bar- code	TCAACC AGCCAG CCT	No
4	12:46321441- 46321441/T/ G-F1	AATCACACTCCA TAGGTATCATTTTC A	gaggtcatacacATGTGA TACTTTGGGAATG AAG	1	Bar- code	CTATGG AGTGTG ATT	No
5	12:46321441- 46321441/T/ G-F2	TTCATTTCATTGT TTAAGATCAGCA	taggacgttcCTTCTGAA CACCAAATTGGAAA	2	Bar- code	AAACAA ATGAAT GAA	No

Template Nucleic Acid

[0103] Samples comprising template nucleic acid molecules to be assayed using the methods disclosed herein can be obtained from a variety of sources including, but not limited to, tissue biopsies, blood draws, buccal swabs, hair, sweat, skin, semen, and mucus. In some embodiments, the sample comprises cells from a subject, for example, circulating tumor cells, blood cells, skin cells, and the like. In some embodiments, the sample comprises cell free nucleic acid, such as, but not limited to, cell free tumor nucleic acid and cell free fetal nucleic acid. In some embodiments, the template nucleic acid molecule is isolated or purified before amplification. Methods of isolating and purifying nucleic acids are well known in the art. Template nucleic acid molecules comprise at least one target nucleic acid sequence. The target sequence is flanked by primer binding sites. In some embodiments, the template is a DNA molecule. In some embodiments, the template is an RNA molecule. In some embodiments, the template may be double-stranded, while in other embodiments, the template is single-stranded.

[0104] In some embodiments, the target nucleic acid is a portion of a gene such as, but not limited to, ABCC8, ABLIM3, ACBD3, ACIN1, ACSL5, ACTA2, ACVR1, ACVR1B, ACVR1C, ACVR2B, ADAMTSL3, ADORA2A, AEBP2, AES, AFAP1, AGAP1, AKR7A2, AKT1, ALK, AMHR2, AMPD3, ANGPTL6, ANO7, APC, APOL2, AQP4-AS1, ARHGEF3, ARID1A, ARIDSA, ARIH1, ARNT, ATM, ATP5A1, ATP9B, ATXN7L1, AX747372, BAG1, BAIAP2L1, BECN2, BMP4, BMP8A, BMP8B, BMPR1A, BMPR1B, C12orf60, C17orf89, C10orf210, C6orf10, C6orf211, C9orf40, CACNA1A, CACNA1H, CACNA2D4, CAMK1D, CAMKMT, CARM1, CAST, CBS, CCBE1, CDC40, CDH23, CDH4, CDKN2B, CHRNA4, CLASP1, CLCA1, CLDN2, CLIC3, CNN3, CNTN1, COL11A2, COL3A1, COL3A2, COL4A1, COL4A5, COL4A6, COL5A1, COL5A2, COL6A2, COL6A3, COX7A2L, CRADD, CREBBP, CRY2, CSGAL-NACT2, CTBP2, CYP2S1, DAG1, DCAF8, DCAF8, DCAF8, DLAT, DLGS, DLGAP4-AS1, DNAH3, DOCK4, DOCK8, DOPEY1, DPYSLS, DYNCH1, DYNCH2, DYRK2, E2F4, E2F6, ECI2, EEF1DP3, EHD4, EIF2B5, EIF4G3, ELAC2, ELK3, EMD, EMX20S, EPPK1, EPT1, ERBB4, ERCCS, ETS2, ETV4, FAM107B, FAM13B, FAM175A, FAM83E, FAV, FBN1, FBN2, FBN3, FBXO28, FGFR2, FHL2, FIRRE, FLNA, FLT3, FOXA3, FOXG1-AS1, FST, GABRG1, GALM, GAPDH, GDF6, GDF7, GLI2, GLI3, GLRXS, GLT8D2, GOLPH3, GPD2, GPR68, GPRASP1, H2AFX, HDAC4, HHAT, HIST1H2AH, HIST2H2AB, HK1, HMCN1, HMSD, HNF4A, HNRNPU, HOXD3, HPS3, HS3ST3A1, IDH1, IFNG, IKBKAP, IMP3, INHBA, INPP4B, INPP5A, IQCK, JAG1, JWT213-1, JWT213-2, JWT213-3, JWT213-4, JWT213-5, JWT213-6, JWT213-7, JWT213-8, JWT213-9, JWT307_1, JWT307_2, JWT307_3, JWT307_4, JWT307_5, JWT307_6, JWT307_7, JWT310-1, JWT310-2, JWT310-3, JWT310-4, JWT310-5, JWT310-6, JWT310-7, JWT311-1, JWT311-2, JWT311-3, JWT311-4, JWT311-5, JWT311-6, JWT311-7, JWT312-1, JWT312-2, JWT312-3, JWT312-4, JWT312-5, JWT312-6, JWT312-7, JWT312-8, JWT312-9, JWT313-1, JWT313-2, JWT313-3, JWT313-4, JWT313-5, JWT313-6, JWT313-7, JWT313-8, JWT313-9, JWT364_1, JWT364_2, JWT364_3, JWT364_4, JWT364_5, JWT364_6, JWT364_7, KANSL1, KCNQ1, KDM3A, KDR, KIRREL3, KLF13,

KLHL14, KMTD2, L3MBTL1, LACTB2, LAMA2, LAMA3, LEFTY1, LINGO4, LMAN2L, LRRC4C, LSAMP, LTBP1, LTBP2, LTBP3, LZTS2, MAD1L1, MAD2L1, MAEA, MAGI2, MAML2, MAP3K7, MAPK1, MAPK3, MAPK8IP2, MARK3, MAT2A, MATR3, MBNL2, MCL1, MCU, MECP2, MED12, MED29, MEF2A, MEGF6, MESD, METTL17, MIER2, MIR181A1HG, MKL1, MKL2, MLH1, MOB2, MPRIP, MRPL32, MRS2, MTCH1, MTOR, MUC16, MUC3A, MYC, MYH11, MYH11, NDE1, MYH11; MYH11, MYLK, MYLK-AS1, MYOCD, NA, NDFIP2, NDUFC1, NEK9, NF1, NFKB1, NGEF, NME4, NME4, DECR2, NOL9, NOTCH1, NOTCH3, NPLOC4, NRG4, NRM, NRTN, NTM, NUCB1, NUDT16, NUDT16L1, OAS3, OR4K3, OSTC, PAG1, PCDH15, PDCD6, PDE4DIP, PDSSA, PHC1, PHF12, PHKG1, PIK3R1, PLEKHG6, PLXDC2, PMM2, POLG2, POLR3B, PPARGC1A, PPHLN1, PPP1R14A, PPP1R15B, PRAF2, PRDM16, PRKG1, PRPH2, PRTG, PTGDR, PTPN12, PTPN14, PTPRC, PTPRS, PUS7, RABL6, RALGAP1, RAPGEF4, RBM10, REPS2, RHBDF2, RIN2, RNF175, RNU1-35P, RNU1-35P, RP11-149P24.1, ROCK1, ROCK2, RPRD2, RSF1, RUSC1, SAFB2, SASH1, SCAF11, SCARF1, SEPT11, SH3GLB2, SHPK, SHPK, SHPK, SHROOM3, SIKE1, SIPA1L1, SIRPA, SK213, SK215, SLAIN1, SLC1A4, SLC25A48, SLC2A10, SLC4A1AP, SLMO2, SLTM, SLX4, SMAD3, SMAD4, SMAD5, SMAD6, SMAD7, SMARCA4, SMLR1, SMTNL1, SMURF1, SNK307, SNK310, SNK311, SNK312, SNK313, SNK364, SNK380, SNK382, SNK383, SNK384, SNK385, SNK386, SOX21-AS1, SOX9, SPOCK2, SPRED1, SPSB2, SRGN, SRP68, SRRM2-AS1, ST6GAL1, STK16, STRN3, SUCLA2, SUCO, SWI5, SYNE2, TAB1, TBC1D13, TBCE, TCERG1, TCF4, TERT, TFB2M, TFDP1, TGFB1, TGFB3, TGFB1, TGFB2, THBS1, TMEFF2, TMEM132C, TMEM2, TMEM268, TNPO1, TPCN2, TPM3, TPRX1, TRAM1, TRAPPC9, TRPM1, TSC2, TSHZ2, TTN, TUBG1, TUBGCP3, TULP4, UBAP2, UBE2I, UBE2W, UHRF1, UNC45A, UNG, UROC1, USP24, USP34, USP8, VANGL1, VIPR2, VPS13D, WDR35, WDR45B, WDR7, WDSUB1, WHSC1, YARS2, YIPF3, ZFH4, ZFYVE16, ZFYVE9, ZMIZ1, ZNF223, ZNF292, ZNF3, ZNF362, ZNF451, ZNF517, ZNF593, ZNF630, ZNRF3, or ZSCAN5A.

[0105] The subject from whom the template nucleic acid molecule sample is obtained can be any organism. In some embodiments, the subject is a vertebrate. In some embodiments, the subject is a mammal such as a human, mouse, rat, dog, cat, horse, cow, sheep, or other domesticated mammal. In some embodiments, the mammal is a human. In some embodiments, the subject from whom the sample is obtained has or is suspected of having a disease or condition associated at least in part with a genetic variant or variants.

Polymerases

[0106] The methods provided herein use a nucleic acid polymerase to amplify a target nucleic acid sequence. Because some polymerases have high error rates (incorporating the wrong nucleotide at a position in a synthesized nucleic acid), selection of a suitable polymerase is an important concern. Sequence errors introduced by a polymerase confound authentic sequence data, making discernment of low frequency variants unreliable or expensive due to the amount of coverage necessary to overcome the polymerase's error rate. High-fidelity polymerases, are par-

ticularly well-suited for use in the presently disclosed methods, and can be used to synthesize copies of a target nucleic acid sequence that potentially harbors a low-frequency variant. Such high-fidelity polymerases introduce fewer nucleotide sequence errors than non-high-fidelity polymerases. Thus, in some embodiments, the nucleic acid amplification reactions comprise a high-fidelity nucleic acid polymerase. For example, in some embodiments, nucleic acid reactions comprise a Phusion high-fidelity DNA polymerase (New England Biolabs (NEB)). This polymerase has a reported error rate of 4.4×10^{-7} errors per base in Phusion HF buffer and 9.5×10^{-7} errors per base in GC buffer. *Thermus aquaticus* (Taq) polymerase has a 50-fold higher error rate than the error rate of the Phusion high-fidelity polymerase. Other polymerases may be used to amplify nucleic acids according to the presently disclosed methods, but an increase in polymerase error rates may decrease the reliability of the method. Table 2 provides a summary of the differences between the high-fidelity Phusion DNA polymerase and the *Pyrococcus furiosus* and the Taq DNA polymerases (HF=high-fidelity; “GC Buffer” refers to a buffer suited for reactions amplifying a target rich in G and/or C). To overcome such errors generated by non-high-fidelity polymerases, additional coverage of the interrogated nucleic acid may be necessary, resulting in increased costs.

TABLE 2

Polymerase Comparison		
Polymerase	1 kb Template	3 kb Template
Phusion High-Fidelity DNA Polymerases (HF Buffer)	1.32%	3.96%
Phusion High-Fidelity DNA Polymerases (GC Buffer)	2.85%	8.55%
<i>Pyrococcus furiosus</i> DNA polymerase	8.4%	25.2%
Taq DNA polymerase	68.4%	>200%

Overview of the Method

[0107] The methods disclosed herein are suitable for detecting low frequency variants. The methods described herein involve detecting the presence or absence of low frequency genetic variation in a nucleic acid molecule by amplifying the nucleic acid sequence of interest using multiple pairs of primers. Each pair of primers comprises a forward primer and a reverse primer, each having a unique binding sequence complementary to a target polynucleotide, wherein the intervening sequences between each pair of primers (i.e., the amplified nucleic acid sequence) at least partially overlap. The resulting overlapping amplicons are sequenced using a Next Generation Sequencing platform, which provides the deep coverage necessary to validate low frequency variants. The sequencing reads are aligned, and determinations regarding the presence or absence of genetic variation are made. The sequencing data can be used for further characterization of any detected genetic variation (i.e., alternative allele fraction).

[0108] In some embodiments, the low frequency variant is a known variant, and the methods disclosed herein may be used to confirm the variant’s presence and/or characteristics (i.e., its alternate allele frequency). In some embodiments, the low frequency variant originated during a germline event, while in other embodiments, the low frequency variant to be interrogated originated during a somatic event.

In some embodiments, the low frequency variant is a silent variant, a missense variant, or a nonsense variant. In some embodiments, the low frequency variant alters a splice site or is an insertion or deletion.

Amplification

[0109] In some embodiments, nucleic acid amplification reactions comprise a template nucleic acid molecule having a target nucleic acid sequence, at least three primer pairs suitable for interrogating the target nucleic acid, nucleotides, and a polymerase. Due to the use of at least three primer pairs in the amplification, the overall method described herein can be referred to a triple-primer PCR sequencing. In some embodiments of the present disclosure, the reaction further comprises a buffer that provides a suitable ionic environment for the polymerase to synthesize a nucleic acid molecule. In some embodiments, the reaction comprises a buffer having essential cofactors (e.g., magnesium) necessary for polymerase function. In some embodiments, the cofactors necessary for proper polymerase function are added to the reaction independently of the buffer.

[0110] In some embodiments, the amplification reaction comprises labeled nucleotides, wherein the labeled nucleotides facilitate efficient capture of any amplicon that comprises one or more labeled nucleotides. Referring to FIG. 3, a nucleotide may be labeled with biotin, and amplicons incorporating the biotin-labeled nucleotides can be captured on streptavidin beads or other media or substrate comprising streptavidin. These captured amplicons can be used as templates for a subsequent amplification reaction, thereby enriching the captured amplicons.

[0111] In some embodiments, separate nucleic acid amplification reactions are prepared for each pair of primers. For example, amplifying a target nucleic acid sequence may comprise at least three reactions according to the methods described herein, wherein each reaction comprises one of three different pairs of primers. The primers, as discussed supra, are used in amplification reactions that generate overlapping amplicons (i.e., semi-redundant interrogation of the target nucleic acid sequence), thereby reducing the probability of impaired detection of variants or skewed downstream determination of alternate allele fractions due to amplification bias. In some embodiments, a single amplification reaction will comprise all pairs of primers. Combining the different primers into a single amplification reaction will generate a greater number of distinct amplicons.

[0112] In some embodiments, the amplification reactions are polymerase chain reactions (PCR). PCR reactions undergo multiple thermocycles, wherein each thermocycle comprises a denaturing step, an annealing step, and an extension step. During the denaturation step, the reaction is incubated at or above 90° C., which is a sufficient temperature, in some embodiments, to cause a double-stranded DNA molecule to denature into single DNA strands or to cause the nucleic acid molecule to undergo a conformational change that is more conducive for an amplification reaction.

[0113] The annealing step comprises complementary binding of the primers to the template nucleic acid and occurs at a lower temperature than that used in the denaturing step. In some embodiments, each primer will be designed to anneal to a complementary nucleic acid sequence at a temperature of between about 50° C. and about 65° C. In some embodiments, the annealing temperature is about 50° C., 51° C., 52° C., 53° C., 54° C., 55° C., 56° C.,

57° C., 58° C., 59° C., 60° C., 61° C., 62° C., 63° C., 64° C., or 65° C. about 50° C., 51° C., 52° C., 53° C., 54° C., 55° C., 56° C., 57° C., 58° C., 59° C., 60° C., 61° C., 62° C., 63° C., 64° C., or 65° C. In some embodiments, the temperature at which the primers anneal to the nucleic acid template can be modified by adjusting conditions (e.g., salt concentration) in the sample or in the amplification reaction. One skilled in the art will understand how changing sample or reaction conditions can affect the temperature at which a primer binds to template nucleic acid.

[0114] In the extension step of a PCR cycle, the primers annealed to the template nucleic acid's primer binding sites are extended by a polymerase to produce a nucleic acid molecule that is complementary to a portion of the template nucleic acid molecule. A proper extension temperature is at or about the optimal temperature for the polymerase to synthesize a nucleic acid molecule. In some embodiments, the extension temperature is between about 65° C. and 75° C. In some embodiments, the extension temperature is about 65° C., 66° C., 67° C., 68° C., 69° C., 70° C., 71° C., 72° C., 73° C., 74° C., or 75° C. In some embodiments, the extension temperature may be 5, 10, 15, 20, or 25% higher or lower than the optimal temperature of the polymerase. Those skilled in the art will understand how to adjust the temperatures, or other reaction conditions, necessary for successful PCR amplification of a nucleic acid sequence.

[0115] In some embodiments, the template nucleic acid is amplified isothermally. For example, helicase dependent amplification is an isothermal amplification method that utilizes a helicase, rather than high temperatures, to separate the strands of a duplex nucleic acid. By not requiring a denaturation step, the isothermal reaction can be incubated at or about the optimal temperature of the polymerase. However, in some embodiments, the isothermal amplification reaction comprises an initial heat denaturation step. Exponential amplification is achieved by incubating the reaction at an isothermal temperature, which obviates the need for thermocycling equipment. Other isothermal amplification techniques are known in the art, and one skilled in the art would understand how to optimize these techniques to comport with the methods described herein.

[0116] Referring to FIGS. 4 and 5, in some embodiments, the amplification reaction products (amplicons) are pooled. This allows simultaneous sequencing of the amplicons generated by the different amplification reactions, which decreases reagent costs and the burden on laboratory personnel and equipment. In some embodiments, the amplification reactions are not pooled prior to sequencing. Pooling, in some embodiments, comprises combining all the amplicons, while in some embodiments, pooling of only a subset of the amplification reactions is required. Additionally, in some embodiments, only a portion of each amplification reaction is pooled, and the remaining unpooled amplification reactions are assayed in parallel with different techniques.

[0117] In some embodiments, the amplification reaction products are purified or isolated before pooling. Methods for isolating and purifying nucleic acids are well known in the art, and there are many commercially available kits for purifying or isolating amplicons. In some embodiments, purifying or isolating amplicons occurs after pooling. In some embodiments, enriched amplicons resulting from biotin:streptavidin capture and reamplification, can be purified using streptavidin to bind and separate all biotin labeled amplicons.

[0118] In some embodiments, the amplicons are assessed prior to being sequenced. Assessing the amplicons can include, for example, gel electrophoresis, real time detection, or spectrophotometric determination of amplicon concentration. For example, amplicons may be assessed using a TapeStation (Agilent) or Bioanalyzer 2100 (Agilent). These analyses allow an investigator to determine if the amplification reaction generated sufficient amounts of high quality amplicons for subsequent sequencing.

Sequencing

[0119] Sequencing of the overlapping amplicons provides multiple independent interrogations of a variant nucleotide or nucleic acid sequence compared to using a single pair of primers. Traditional Sanger sequencing platforms can be used to sequence the overlapping amplicons, but this approach is inefficient for detecting rare variants. Conversely, Next Generation Sequencing (NGS) platforms can generally accommodate thousands of sequencing reactions run in parallel, thereby providing deeper coverage than is possible with Sanger sequencing. For example, referring to FIG. 6, the Ion Torrent system can generate nearly twenty million reads with 93% ion sphere particle (ISP) loading. Ion sphere particles used in the Ion Torrent system are conjugated directly or indirectly to a nucleic acid comprising the sequence of interest adjacent to a nucleic acid sequence complementary to the adapter described supra. In detecting, characterizing, or validating low frequency variants, this increased coverage enables distinguishing true variants from errors introduced during amplification, sequencing, or data processing.

[0120] The amplicons to be sequenced are, by design, generally less than 300 nucleotides in length, and there are several NGS platforms that can cost-effectively generate sequencing data at the desired coverage level. For example, ThermoFisher's Ion Torrent and Illumina's MiSeq can each generate maximum read lengths of approximately 250 nucleotides. Other NGS approaches are available for shorter or longer read lengths. For example, Illumina's HiSeq platform has a maximum read length of about 150 nucleotides, while the Roche 454 platform can generate at least 400 nucleotide reads. One skilled in the art will be to determine which platform can be used to generate the desired sequencing data, and will optimize the adapters on each primer to comport with that platform.

Data Processing and Analysis

[0121] In some embodiments, the sequencing data is assessed for quality before alignment, and those reads not possessing the required quality characteristics are removed from the data set. Typically, quality control of sequencing reactions comprises establishing a signal-to-noise threshold, and reads that do not meet the threshold are discarded. Such quality control lessens the probability of erroneous base calls in a read that would decrease reliability of the assay.

[0122] Sequencing data generated using the disclosed methods can be processed to accurately determine alternate allele frequencies. Referring to FIG. 7, in some embodiments, the sequencing data is first demultiplexed by grouping together all reads having the same index sequence. Each pair of primers used to amplify a target nucleic acid sequence has a unique index sequence, such that data generated for the products of distinct amplification reactions

will be segregated into distinct bins based on their index sequence. All sequences having the same index sequence will be binned together and segregated from sequences having different index sequences. This demultiplexing of the sequencing data allows for three independent determinations of the alternate allele fraction for variants detected in the target nucleic acid sequence and the assignments of confidence intervals. In some embodiments, the average alternate allele fraction is determined by averaging the three individual alternate allele fractions.

[0123] The data in each bin is aligned to provide maximal sequence identity between the individual reads. For example, if a read has a single nucleotide deletion, the alignment will incorporate the deletion into the read's aligned sequence so that the nucleotide sequences on either side of the deletion align with other reads that do not have the deletion. Referring to FIG. 8, indels are elevated in Ion Torrent sequencing, and these errors can mask true alleles (especially low frequency variants) (top panel). However, the Pullox Algorithm can identify and correct about 97% of such indel errors and does not impact mosaic alleles (middle panel). This program can also reduce background noise up to 50%. The processed data can be mapped to the genome or template nucleic acid and is able to identify the target allele (bottom panel).

[0124] Primer binding sites are also identified (FIG. 9) and removed from the sequencing data. Because these sequences are known, they can be readily identified and removed, which avoids analyzing possible false positive and false negative results in these sequences.

[0125] In some embodiments, all but one read having the same unique molecular identifiers will be removed from the data set, which indicates multiple amplification reactions that used the exact same primer. These duplicated amplifications reactions are not considered independent interrogations of the nucleic acid. Retention of such redundant data could impact alternate allele fraction determination. In some embodiments, accurate determination or validation of alternate allele frequencies of about 0.025% comprise removing redundant reads from the data. In some embodiments, wherein the alternate allele fraction is known to be 0.1% or greater, removal of redundant reads may not be necessary due to the deep coverage available in Next Generation Sequencing platforms. Once the alignment is set in each bin, the alternate allele frequencies for variants in each bin are determined.

[0126] The methods provided can distinguish between germline and somatic events resulting in genetic variation. Referring to FIG. 10A, a genetic variant derived from a germline event, which should approach an alternative allele frequency of about 50%, is shown. Three panels of sequencing data are separated by the large shaded boxes, wherein each panel presents a subset of sequencing data for amplicons generated from different amplification reactions. In each panel, the allele frequency is nearly identical in each panel (Panel 1: 49.5% (112,000× coverage); Panel 2: 49.9% (75,000× coverage); and Panel 3: 50.0% (126,000× coverage). The alternate allele frequencies are then averaged for each variant and a confidence interval assigned. Those skilled in the art will understand how the frequencies are determined and will know that commercially available algorithms can be employed.

[0127] A somatic event occurring in a single subject will likely have a much lower allele frequency than an inherited

allele, and a subject having a genetic variant derived from a somatic event is said to be mosaic for the variant. As shown in Table 3, the alternate allele frequencies (AAF) observed in three different amplicon samples are about 1%, well below the frequency expected in an individual for an inherited allele, which suggests the variant is a somatic mosaic variant. For example, for the sequencing reads of amplicons generated using the Primer 1 set of primers, 416 reads out of 37,779 total reads contained the alternate allele (FIG. 10B). The “Background AAF” is the alternative allele frequency of variants detected in the regions flanking the alternate allele (also referred to as the “background rate”). In some embodiments, sequencing data of the primer binding sites is removed prior to determining a background rate. This improves the accuracy of the background rate because sequencing errors are more prevalent for regions near the adapter binding sites (e.g., primer binding sites).

TABLE 3

Alternative Allele Fractions			
Primer #:	Allele Counts	AAF	Background AAF
Primer 1	416/37779	1.09%	0.0009%
Primer 2	123/13064	0.94%	0.0045%
Primer 3	529/50141	1.04%	0.0027%
Average	—	1.02% ± 0.19% (p = 0.0009)	0.0025%

Method Comparison

[0128] Two methods are currently used to detect and quantify rare variants, droplet digital PCR (ddPCR) and Sanger sequencing of TOPO (Topoisomerase-based) cloned nucleic acids. Referring to Table 4, the estimated cost of the method described herein (“mosaic validation method”) is about 90% less expensive than ddPCR and 85× less expensive than the Sanger sequencing/TOPO cloning method. Furthermore, the Sanger sequencing/TOPO cloning method is much less sensitive as its lowest level of reliable detection is an alternate allele fraction of 0.5%. While the purported resolution of ddPCR is an alternate allele fraction of 0.1%, it is not reliable for alternate allele fractions of 0.02% that are within the reliable range of the presently disclosed methods.

[0129] Additionally, high-throughput Next Generation Sequencing platforms used in the presently disclosed methods can run massive parallel reactions. Conversely, both Sanger Sequencing/TOPO cloning and ddPCR have relatively limited throughput, thereby increasing cost and time requirements. ddPCR, while having higher throughput than the Sanger sequencing/TOPO cloning method, does not enjoy the throughput of the presently described methods. Additionally, ddPCR primers are labeled with a relatively expensive fluorophore.

TABLE 4

Method Comparison			
	ddPCR	Mosaic Validation Method	Sanger + TOPO Cloning
Estimated Cost to Validate allele	\$256	\$35	\$3,004

TABLE 4-continued

	Method Comparison		
	ddPCR	Mosaic Validation Method	Sanger + TOPO Cloning
Cost of Amplification Primers	\$250 (1 set)	\$27 (3 sets)	\$4 (1 set)
Cost of Sequencing/Amplification Resolution Throughput	\$6/triplicate 0.1% AAF Low-medium	\$8/3 primers 0.02% AAF High	\$3,000/mutation (1,000 colonies at \$3 per colony) 0.5% AAF Low

Detecting and Monitoring Disease

[0130] The methods described herein can be used for the detection and/or monitoring of a disease. The detection and characterization of disease-associated variants, including somatic mosaic variants, can provide information relevant for diagnosing a disease, determining the progression or regression of disease, and treating disease. For example, when a cancer cell arises after a somatic event, or when circulating tumor cells are present in a subject, the methods described herein can be used to detect of these cells.

[0131] A subject having a disease may undergo periodic testing to determine if the number of a diseased cells is increasing, decreasing, or static. For example, a subject that has cancer may determine the alternative allele frequency of a cancer marker present in samples after the cancer is detected or after treatment has begun. Changes in the alternative allele frequency of the cancer marker would indicate a change in the number of cells carrying the marker (e.g., cancer cells) present in the sample. If the alternative allele frequency is greater than that observed in a previous sample, the subject's cancer is likely progressing or not responding effectively to treatment. If the alternative allele frequency remains static relative to an earlier sample, the disease may be responding treatment sufficiently to stop disease progression, but perhaps not to a level sufficient for disease regression or remission. If the alternative allele frequency decreases relative to an earlier sample, the subject's disease may be regressing, and the absence of such cells (i.e., AAF=0) may signify remission.

Kits and Compositions for Detecting and Characterizing Low Frequency Genetic Variation

[0132] In another embodiment, kits and compositions are provided that advantageously allow for the detection and/or quantification of the presence of low frequency genetic variation in a subject sample (e.g., blood or serum). In one embodiment, the kit includes a composition comprising reagents for performing an amplification reaction, including multiple pairs of forward and reverse primers as described herein. In some embodiments, the reagents include nucleotides, labeled nucleotides, a buffer, a cofactor, and/or a polymerase. In some embodiments, the kit comprises a sterile container that contains the amplification reaction reagents; such containers can be boxes, ampoules, bottles, vials, tubes, bags, pouches, blister-packs, or other suitable container forms known in the art. Such containers can be made of plastic, glass, laminated paper, metal foil, or other materials suitable for holding amplification reagents.

[0133] In one embodiment, the kit comprises high-quality (PAGE-purified) RNA or DNA-based primers, premixed at proper concentrations. In some embodiments, the kit comprises reagents for biotin labeling for higher sensitivity assays. In some embodiments, the kit comprises a preselected polymerase (e.g., Phusion U if using RNA primers, or another option) with high fidelity (100x improved error rates compared to a reference polymerase (Taq polymerase). In some embodiments, the kit comprises duplicate primers with differing barcodes for testing case/control samples side-by-side. In some embodiments, the kit comprises preselected primers to avoid other mutation sites, non-overlapping binding sites, and the like. In some embodiments, the kit comprises control DNA (e.g., for negative controls). In some embodiments, the kit comprises ddPCR probes for performing ddPCR and sequencing from the same reaction—(i.e., to obtain copy/expression values and genotype correlation).

[0134] In another embodiment, the kit includes a composition comprising reagents for performing a sequencing reaction, including nucleic acid molecules that can specifically bind to an adapter as described above. The reagents, in some embodiments, include nucleotides, labeled nucleotides, a buffer, a cofactor, ion spheres comprising the nucleic acid molecule to be sequenced, and/or enzymes for catalyzing the sequencing reaction. In some embodiments, the kit comprises a sterile container that contains the sequencing reaction reagents; such containers are described above.

[0135] In some embodiments, the kit comprises compositions for amplification and sequencing as described above. Kits may also include instructions for performing the reactions.

[0136] The practice of the present disclosure teaches, unless otherwise indicated, conventional techniques of molecular biology (including recombinant techniques), microbiology, cell biology, biochemistry and immunology, which are well within the purview of the skilled artisan. Such techniques are explained fully in the literature, such as, "Molecular Cloning: A Laboratory Manual", second edition (Sambrook, 1989); "Oligonucleotide Synthesis" (Gait, 1984); "Animal Cell Culture" (Freshney, 1987); "Methods in Enzymology" "Handbook of Experimental Immunology" (Weir, 1996); "Gene Transfer Vectors for Mammalian Cells" (Miller and Calos, 1987); "Current Protocols in Molecular Biology" (Ausubel, 1987); "PCR: The Polymerase Chain Reaction", (Mullis, 1994); "Current Protocols in Immunology" (Coligan, 1991). These techniques are applicable to the production of the polynucleotides and polypeptides of the disclosure, and, as such, may be considered in making and practicing the compositions and methods disclosed herein. Particularly useful techniques for particular embodiments will be discussed in the sections that follow.

[0137] The following examples are put forth to provide those of ordinary skill in the art with a complete disclosure and description of how to perform the amplification, sequencing, and quantifying methods presently disclosed, and are not intended to limit the scope of what the inventors regard as their invention.

EXAMPLES

Example 1: Detecting Alleles with an Alternate Allele Fraction (AAF) at or Above 0.1%

[0138] To identify low frequency genetic variation in a target nucleic acid sequence with an alternate allele fraction

(AAF) of 0.1% or greater, three pairs of primers were designed to yield overlapping amplicons. Each pair of primers comprised a forward and a reverse primer, with each primer having a nucleotide sequence complementary to a portion of the target nucleic acid sequence. Each primer had an adapter at or near its 5' terminus and upstream from its complementary nucleic acid sequence. The adapter's nucleic acid sequence was complementary to a nucleic acid sequence used in a Next Generation Sequencing (NGS) platform, such as Ion Torrent or Illumina's MiSeq. Additionally, the reverse primer for each pair of primers further comprised an index sequence upstream from the primer's complementary nucleic acid sequence that was unique to the pair.

[0139] Three distinct amplification reactions were prepared, each comprising one of the three pairs of primers. The reactions comprised 1.0 μ M primers, 1 $\times$ final concentration of 5 $\times$ Phusion High-fidelity Buffer (NEB), 200 μ M dNTPs, 1.0 units of Phusion High-fidelity Polymerase (NEB), and about 25 to 50 ng of template DNA. The reactions were subjected to an initial denaturation step of 30 seconds at 98° C. followed by 20 cycles of 98° C. (denaturing the template DNA) for 10 seconds, 62° C. (annealing the primers to the template nucleic acid) for 20 seconds, and 72° C. (to extend the DNA product) for 30 seconds. After cycling, the reactions were subjected to an additional 10 minutes at 72° C. as a final extension step.

[0140] 5 μ l of each PCR product were then pooled and purified using a ThermoFisher MagJet purification kit (any kit that removes products <100 base pairs in length can be used). The purified reaction products were resuspended in 20 μ l of water, mixed, and incubated for two minutes. The reactions were then placed on a magnet for two minutes, and the eluted DNA was removed. About 1 μ l was run on a TapeStation or a Bioanalyzer 2100 to confirm quality.

[0141] Aliquots of the amplicons generated from a single round of amplification were analyzed on a Bioanalyzer 2100. This amplification strategy yielded detectable amplicons at the expected time point (i.e., between 50 and 60 seconds for the control (FIGS. 11A and 11B) and between 70 and 80 seconds for the amplification performed according to the single round amplification methods described herein (FIGS. 11A and C)). The dark bands at approximately 43 and 113 seconds are control nucleic acids. PicoGreen (ThermoFisher) is then used to measure the concentration of the PCR product, which was subsequently diluted to 100 pM.

[0142] The purified PCR reaction products were sequenced using the Ion Torrent system (ThermoFisher Scientific) to generate sequencing reads that comprise the nucleic acid sequence of the target nucleic acid. The sequencing reads were demultiplexed, or segregated, into different bins depending on the detected index sequence. Table 5 provides a summary of the observed alternate allele fractions detected using this method.

TABLE 5

Observed alternate allele fractions

PrimerID	Chr	AlleleStart	Ref	Alt	Gene	IT Read Depth	Alt Allele Depth	Background AAF (within 50 ms)	Stdev AAF	Variance Background AAF	Stdev of Average Background	Confidence interval of Background	Average AAF
2FLNA_X_153579448_A_G_P44201_2	X	153579448	A	G	FLNA	52876	331	1.02173E-05	2.84343E-05	7.9943E-10	6.34682E-05	0.000157664	0.006458143
4SCAF11_12_46321441_T_G_P44201_1	12	46321441	T	G	SCAF11	37184	129	5.45416E-06	1.50886E-05	2.2511E-10	0.00011421	0.000283714	0.002795873
5SCAF11_12_46321441_T_G_P44201_2	12	46321441	T	G	SCAF11	30037	49	3.91357E-05	0.00010047	9.9614E-09	0.00011421	0.000283714	0.002795873
6SCAF11_12_46321441_T_G_P44201_3	12	46321441	T	G	SCAF11	64191	211	4.41917E-05	0.000171234	2.8945E-08	0.00011421	0.000283714	0.002795873
10SLX4_16_3639306_G_A_P44201_1	16	3639306	G	A	SLX4	45836	265	2.45568E-05	5.14068E-05	2.6139E-09	5.74749E-05	0.000142776	0.003440733
11SLX4_16_3639306_G_A_P44201_2	16	3639306	G	A	SLX4	52791	145	4.46418E-05	7.51276E-05	5.5638E-09	5.74749E-05	0.000142776	0.003440733
12SLX4_16_3639306_G_A_P44201_3	16	3639306	G	A	SLX4	41805	75	1.69855E-05	4.18365E-05	1.7304E-09	5.74749E-05	0.000142776	0.003440733
16LAMA3_18_21453038_C_T_P44201_1	18	21453038	C	T	LAMA3	44807	167	1.75378E-05	5.46427E-05	2.9554E-09	8.21671E-05	0.000204114	0.0038759
17LAMA3_18_21453038_C_T_P44201_2	18	21453038	C	T	LAMA3	10076	37	3.19472E-05	9.90798E-05	9.6703E-09	8.21671E-05	0.000204114	0.0038759
18LAMA3_18_21453038_C_T_P44201_3	18	21453038	C	T	LAMA3	46352	196	6.83073E-05	8.78932E-05	7.6286E-09	8.21671E-05	0.000204114	0.0038759
19FLNA_X_153587777_G_C_P44201_1	X	153587777	G	C	FLNA	48596	289	2.0029E-05	4.94376E-05	2.4189E-09	7.59331E-05	0.000188628	0.0055221
20FLNA_X_153587777_G_C_P44201_2	X	153587777	G	C	FLNA	51421	304	4.46736E-05	0.000116647	1.3412E-08	7.59331E-05	0.000188628	0.0055221
21FLNA_X_153587777_G_C_P44201_3	X	153587777	G	C	FLNA	35689	168	1.85615E-05	3.84922E-05	1.4664E-09	7.59331E-05	0.000188628	0.0055221
28SCAF11_12_46321441_T_G_P44201_1	12	46321441	T	G	SCAF11	72595	141	5.63756E-06	1.88736E-05	3.521E-10	6.8039E-05	0.000169018	0.001519564
29SCAF11_12_46321441_T_G_P44201_2	12	46321441	T	G	SCAF11	17298	17	3.91321E-05	0.000105292	1.0939E-08	6.8039E-05	0.000169018	0.001519564
30SCAF11_12_46321441_T_G_P44201_3	12	46321441	T	G	SCAF11	60601	99	1.62602E-05	5.56879E-05	3.0695E-09	6.8039E-05	0.000169018	0.001519564
34SLX4_16_3639306_G_A_P44201_1	16	3639306	G	A	SLX4	71852	100	2.80354E-05	0.000109029	1.1752E-08	9.16386E-05	0.000227643	0.001528917
35SLX4_16_3639306_G_A_P44201_2	16	3639306	G	A	SLX4	20195	43	2.27701E-05	0.000109229	1.1755E-08	9.16386E-05	0.000227643	0.001528917
36SLX4_16_3639306_G_A_P44201_3	16	3639306	G	A	SLX4	44100	47	1.88841E-05	4.12854E-05	1.6851E-09	9.16386E-05	0.000227643	0.001528917
40LAMA3_18_21453038_C_T_P44201_1	18	21453038	C	T	LAMA3	119239	348	1.75474E-05	5.56879E-05	3.0695E-09	6.79599E-05	0.000168822	0.001625798
41LAMA3_18_21453038_C_T_P44201_2	18	21453038	C	T	LAMA3	44385	27	2.89431E-05	7.37946E-05	5.3721E-09	6.79599E-05	0.000168822	0.001625798
42LAMA3_18_21453038_C_T_P44201_3	18	21453038	C	T	LAMA3	89592	121	3.1764E-05	7.40446E-05	5.4141E-09	6.79599E-05	0.000168822	0.001625798
43FLNA_X_153587777_G_C_P44201_1	X	153587777	G	C	FLNA	53971	238	2.50405E-05	6.37295E-05	4.0196E-09	7.02955E-05	0.000174624	0.003419407
44FLNA_X_153587777_G_C_P44201_2	X	153587777	G	C	FLNA	70189	280	3.09214E-05	8.27202E-05	6.7489E-09	7.02955E-05	0.000174624	0.003419407
45FLNA_X_153587777_G_C_P44201_3	X	153587777	G	C	FLNA	35499	66	2.36192E-05	6.40637E-05	4.0559E-09	7.02955E-05	0.000174624	0.003419407
46ZNF223_19_44571260_C_A_P44201_1	19	44571260	C	A	ZNF223	26856	0	1.20017E-05	2.67638E-05	7.0851E-10	9.09222E-05	0.000225863	0.000220546
47ZNF223_19_44571260_C_A_P44201_2	19	44571260	C	A	ZNF223	36859	0	1.57225E-05	7.08468E-05	4.9557E-09	9.09222E-05	0.000225863	0.000220546
48ZNF223_19_44571260_C_A_P44201_3	19	44571260	C	A	ZNF223	37785	25	6.62844E-05	0.000139445	1.9136E-08	9.09222E-05	0.000225863	0.000220546
49FLNA_X_153579448_A_G_P44201_1	X	153579448	A	G	FLNA	50890	95	9.53088E-06	2.75309E-05	7.5006E-10	6.09297E-05	0.000151358	0.002228303
50FLNA_X_153579448_A_G_P44201_2	X	153579448	A	G	FLNA	22626	61	1.76593E-05	3.42833E-05	1.1615E-09	6.09297E-05	0.000151358	0.002228303
58SLX4_16_3639306_G_A_P44201_1	16	3639306	G	A	SLX4	56582	42	2.07782E-05	7.67312E-05	5.82E-09	7.90631E-05	0.000196404	0.001099338
59SLX4_16_3639306_G_A_P44201_2	16	3639306	G	A	SLX4	20626	11	2.99415E-05	0.000102898	1.0412E-08	7.90631E-05	0.000196404	0.001099338
60SLX4_16_3639306_G_A_P44201_3	16	3639306	G	A	SLX4	17306	35	2.07858E-05	5.05071E-05	2.5213E-09	7.90631E-05	0.000196404	0.001099338
67FLNA_X_153587777_G_C_P44201_1	X	153587777	G	C	FLNA	104074	329	1.52718E-05	3.30872E-05	1.0835E-09	6.50076E-05	0.000161488	0.00207533
68FLNA_X_153587777_G_C_P44201_2	X	153587777	G	C	FLNA	30969	60	6.13577E-05	9.73719E-05	9.3284E-09	6.50076E-05	0.000161488	0.00207533
69FLNA_X_153587777_G_C_P44201_3	X	153587777	G	C	FLNA	64753	73	2.66973E-05	4.78617E-05	2.2661E-09	6.50076E-05	0.000161488	0.00207533
70ZNF223_19_44571260_C_A_P44201_1	19	44571260	C	A	ZNF223	21369	0	1.01326E-05	2.94338E-05	8.5693E-10	5.64567E-05	0.000140246	2.34169E-05
71ZNF223_19_44571260_C_A_P44201_2	19	44571260	C	A	ZNF223	22286	1	1.64077E-05	5.01574E-05	2.4839E-09	5.64567E-05	0.000140246	2.34169E-05
72ZNF223_19_44571260_C_A_P44201_3	19	44571260	C	A	ZNF223	39402	1	2.85368E-05	7.94613E-05	6.2212E-09	5.64567E-05	0.000140246	2.34169E-05
82SLX4_16_3639306_G_A_P44201_1	16	3639306	G	A	SLX4	38589	18	2.49275E-05	5.3876E-05	2.8408E-09	6.49472E-05	0.000161338	0.004000115
83SLX4_16_3639306_G_A_P44201_2	16	3639306	G	A	SLX4	46474	20	3.26162E-05	8.51179E-05	7.1416E-09	6.49472E-05	0.000161338	0.004000115
84SLX4_16_3639306_G_A_P44201_3	16	3639306	G	A	SLX4	46122	14	2.9696E-05	5.19921E-05	2.6721E-09	6.49472E-05	0.000161338	0.004000115
91FLNA_X_153587777_G_C_P44201_1	X	153587777	G	C	FLNA	57104	2	1.93602E-05	4.88796E-05	2.3646E-09	5.3862E-05	0.000133801	0.00083567
92FLNA_X_153587777_G_C_P44201_2	X	153587777	G	C	FLNA	77610	5	5.43539E-05	5.72273E-05	3.2301E-09	5.3862E-05	0.000133801	0.00083567
93FLNA_X_153587777_G_C_P44201_3	X	153587777	G	C	FLNA	33644	81	2.4105E-05	5.60543E-05	3.1087E-09	5.3862E-05	0.000133801	0.00083567

TABLE 5-continued

Observed alternate allele fractions														
PrimerID	Chr	AlleleStart	Ref	Alt	Gene	IT Read Depth	Background (within 50 nts)			Stdev Background AAF	Variance Background AAF	Stdev of Average Background	Confidence interval of Background	Average AAF
							Alt Allele AAF	Depth	Stdev					
94ZNF223_19_44571260_C.A.PH4201_1	19	44571260	C	A	ZNF223	50984	1	1.44169E-05	2.46323E-05	6.0016E-10	5.68753E-05	0.000141286	0.000110502	
95ZNF223_19_44571260_C.A.PH4201_2	19	44571260	C	A	ZNF223	85847	24	5.09988E-05	8.53336E-05	7.1731E-09	5.68753E-05	0.000141286	0.000110502	
96ZNF223_19_44571260_C.A.PH4201_3	19	44571260	C	A	ZNF223	30935	1	2.62988E-05	4.43728E-05	1.9311E-09	5.68753E-05	0.000141286	0.000110502	
97FLNA_X_153579448_A.G.PH4201_1	X	153579448	A	G	FLNA	62892	1	1.12555E-05	3.17585E-05	9.982E-10	4.03517E-05	0.000100239	2.46952E-05	
98FLNA_X_153579448_A.G.PH4201_2	X	153579448	A	G	FLNA	68746	4	1.07011E-05	2.88479E-05	8.2285E-10	4.03517E-05	0.000100239	2.46952E-05	
99FLNA_X_153579448_A.G.PH4201_3	X	153579448	A	G	FLNA	27140	0	2.40557E-05	5.6478E-05	3.0637E-09	4.03517E-05	0.000100239	2.46952E-05	
100SCAF11_12_46321441_T.G.PH4201_1	12	46321441	T	G	SCAF11	65401	11	5.89969E-06	2.71282E-05	7.2767E-10	0.000103346	0.000256727	6.27575E-05	
101SCAF11_12_46321441_T.G.PH4201_2	12	46321441	T	G	SCAF11	24696	0	3.15284E-05	8.72093E-05	7.5067E-09	0.000103346	0.000256727	6.27575E-05	
102SCAF11_12_46321441_T.G.PH4201_3	12	46321441	T	G	SCAF11	49802	1	4.07656E-05	0.000155257	2.3807E-08	0.000103346	0.000256727	6.27575E-05	
106SLX4_16_3639306_G.A.PH4201_1	16	3639306	G	A	SLX4	60556	41	2.11937E-05	6.02461E-05	3.5901E-09	5.85276E-05	0.000145391	0.000616922	
107SLX4_16_3639306_G.A.PH4201_2	16	3639306	G	A	SLX4	85121	37	2.53988E-05	7.16005E-05	5.0617E-09	5.85276E-05	0.000145391	0.000616922	
108SLX4_16_3639306_G.A.PH4201_3	16	3639306	G	A	SLX4	33828	25	1.85644E-05	4.05368E-05	1.6246E-09	5.85276E-05	0.000145391	0.000616922	
112LAMA3_18_21453038_C.T.PH4201_1	18	21453038	C	T	LAMA3	141247	17	1.81195E-05	4.55462E-05	2.0533E-09	0.000187047	0.000464651	0.000155433	
114LAMA3_18_21453038_C.T.PH4201_3	18	21453038	C	T	LAMA3	106954	37	2.62147E-05	6.57093E-05	4.2637E-09	0.000187047	0.000464651	0.000155433	
115FLNA_X_153587777_G.C.PH4201_1	X	153587777	G	C	FLNA	48712	0	1.33842E-05	4.73444E-05	2.2184E-09	6.18101E-05	0.000153545	0.00046184	
116FLNA_X_153587777_G.C.PH4201_2	X	153587777	G	C	FLNA	14435	2	2.92084E-05	6.26768E-05	3.8746E-09	6.18101E-05	0.000153545	0.00046184	
117FLNA_X_153587777_G.C.PH4201_3	X	153587777	G	C	FLNA	34613	0	2.62791E-05	7.36629E-05	5.3685E-09	6.18101E-05	0.000153545	0.00046184	
118ZNF223_19_44571260_C.A.PH4201_1	19	44571260	C	A	ZNF223	50603	1	1.31139E-05	3.22556E-05	1.0297E-09	7.24504E-05	0.000179977	6.58723E-06	
119ZNF223_19_44571260_C.A.PH4201_2	19	44571260	C	A	ZNF223	42129	0	2.51869E-05	0.000102637	1.0399E-08	7.24504E-05	0.000179977	6.58723E-06	
120ZNF223_19_44571260_C.A.PH4201_3	19	44571260	C	A	ZNF223	76059	0	3.28129E-05	6.6232E-05	4.3181E-09	7.24504E-05	0.000179977	6.58723E-06	
124SCAF11_12_46321441_T.G.PH4201_1	12	46321441	T	G	SCAF11	81594	0	7.56131E-06	2.77648E-05	7.6222E-10	6.80237E-05	0.00016898	4.57173E-06	
125SCAF11_12_46321441_T.G.PH4201_2	12	46321441	T	G	SCAF11	82193	0	3.60845E-05	8.55122E-05	7.2174E-09	6.80237E-05	0.00016898	4.57173E-06	
126SCAF11_12_46321441_T.G.PH4201_3	12	46321441	T	G	SCAF11	72912	1	3.25164E-05	7.73035E-05	5.9021E-09	6.80237E-05	0.00016898	4.57173E-06	
130SLX4_16_3639306_G.A.PH4201_1	16	3639306	G	A	SLX4	27339	10	3.40321E-05	7.19695E-05	5.1227E-09	6.33819E-05	0.000157449	0.000283253	
131SLX4_16_3639306_G.A.PH4201_2	16	3639306	G	A	SLX4	67412	31	2.06402E-05	6.86263E-05	4.6484E-09	6.33819E-05	0.000157449	0.000283253	
132SLX4_16_3639306_G.A.PH4201_3	16	3639306	G	A	SLX4	41457	1	2.35701E-05	4.80335E-05	2.2807E-09	6.33819E-05	0.000157449	0.000283253	
136LAMA3_18_21453038_C.T.PH4201_1	18	21453038	C	T	LAMA3	116508	90	1.83925E-05	5.41278E-05	2.8999E-09	5.91213E-05	0.000146865	0.00031725	
137LAMA3_18_21453038_C.T.PH4201_2	18	21453038	C	T	LAMA3	91222	4	2.2096E-05	6.01332E-05	3.5735E-09	5.91213E-05	0.000146865	0.00031725	
138LAMA3_18_21453038_C.T.PH4201_3	18	21453038	C	T	LAMA3	59074	8	2.90008E-05	6.37446E-05	4.0126E-09	5.91213E-05	0.000146865	0.00031725	
139FLNA_X_153587777_G.C.PH4201_1	X	153587777	G	C	FLNA	48033	3	2.29118E-05	5.18672E-05	2.6625E-09	5.62488E-05	0.00013973	3.21058E-05	
140FLNA_X_153587777_G.C.PH4201_2	X	153587777	G	C	FLNA	61361	0	3.44416E-05	7.64576E-05	5.761E-09	5.62488E-05	0.00013973	3.21058E-05	
141FLNA_X_153587777_G.C.PH4201_3	X	153587777	G	C	FLNA	29533	1	1.53183E-05	3.28539E-05	1.0683E-09	5.62488E-05	0.00013973	3.21058E-05	
190NA_5_173266954_G.A.A-pancreas_1	5	173266954	G	A	NA	106354	134	6.35041E-05	0.000307	9.3213E-08	0.000187713	0.000466304	0.000853995	
196NA_5_173266954_G.A.A-pancreas_2	5	173266954	G	A	NA	122112	2681	3.94563E-05	9.45171E-05	8.8128E-09	0.000117275	0.000291328	0.020158133	
199NA_5_173266954_G.A.A-pancreas_3	5	173266954	G	A	NA	93898	51	4.59757E-05	0.00010127	1.0122E-08	0.000187713	0.000466304	0.000853995	
205NA_5_173266954_G.A.A-pancreas_4	5	173266954	G	A	NA	39129	799	6.47892E-05	0.000119659	1.3993E-08	0.000117275	0.000291328	0.020158133	
208NA_5_173266954_G.A.A-pancreas_5	5	173266954	G	A	NA	51390	39	3.69912E-05	4.92248E-05	2.3726E-09	0.000187713	0.000466304	0.000853995	
212NA_11_49854989_C.T.A-17_3	11	49854989	C	T	NA	24985	4	2.43233E-05	6.79798E-05	4.5716E-09	5.29008E-05	0.000131413	0.000800048	
213NA_11_49854989_C.T.A-17_1	11	49854989	C	T	NA	95864	0	1.34233E-05	3.22146E-05	1.0254E-09	5.29008E-05	0.000131413	0.000800048	
5SCAF11_12_46321441_T.G.PH4201_1	12	46321441	T	G	SCAF11	42588	20869	3.69888E-06	1.05458E-05	1.0996E-10	9.31199E-05	0.000231323	0.491051667	
5SCAF11_12_46321441_T.G.PH4201_2	12	46321441	T	G	SCAF11	13886	6743	6.22253E-05	0.000150646	2.2396E-08	9.31199E-05	0.000231323	0.491051667	
5SCAF11_12_46321441_T.G.PH4201_3	12	46321441	T	G	SCAF11	54414	27073	2.36841E-05	5.95921E-05	3.5084E-09	9.31199E-05	0.000231323	0.491051667	
10SLX4_16_3639306_G.A.PH4201_1	16	3639306	G	A	SLX4	32196	18299	2.91067E-05	6.47837E-05	4.1528E-09	6.1556E-05	0.000152914	0.542452333	
11SLX4_16_3639306_G.A.PH4201_2	16	3639306	G	A	SLX4	22986	12241	2.16324E-05	6.40091E-05	4.0453E-09	6.1556E-05	0.000152914	0.542452333	
12SLX4_16_3639306_G.A.PH4201_3	16	3639306	G	A	SLX4	47406	24957	2.72861E-05	5.66025E-05	3.1694E-09	6.1556E-05	0.000152914	0.542452333	
19FLNA_X_153587777_G.C.PH4201_1	X	153587777	G	C	FLNA	44686	20278	2.41014E-05	5.75831E-05	3.282E-09	9.8689E-05	0.000245157	0.485348667	
20FLNA_X_153587777_G.C.PH4201_2	X	153587777	G	C	FLNA	64553	33761	5.76983E-05	0.000153505	2.3241E-08	9.8689E-05	0.000245157	0.485348667	

TABLE 5-continued

Observed alternate allele fractions													
PrimerID	Chr	AlleleStart	Ref	Alt	Gene	IT Read Depth	Alt Allele Depth	Background AAF (within 50 nts)	Stdev AAF	Variance Background AAF	Stdev of Average Background	Confidence interval of Background	Average AAF
21FLNA_X_153587777_G_C.PH4201_3	X	153587777	G	C	FLNA	38524	18463	2.10532E-05	5.21856E-05	2.6956E-09	9.8689E-05	0.000245157	0.485348667
25FLNA_X_153579448_A_G.PH4201_1	X	153579448	A	G	FLNA	114923	33468	9.44768E-06	2.42237E-05	5.8074E-10	2.96956E-05	7.3768E-05	0.213852
26FLNA_X_153579448_A_G.PH4201_2	X	153579448	A	G	FLNA	91936	21769	1.45941E-05	3.33324E-05	1.0994E-09	2.96956E-05	7.3768E-05	0.213852
27FLNA_X_153579448_A_G.PH4201_3	X	153579448	A	G	FLNA	38714	4396	1.18446E-05	3.1243E-05	9.654E-10	2.96956E-05	7.3768E-05	0.213852
28SCAF11_12_46321441_T_G.PH4201_1	12	46321441	T	G	SCAF11	89492	11917	4.87432E-05	8.54665E-05	7.2225E-09	8.6509E-05	0.0002149	0.132397333
29SCAF11_12_46321441_T_G.PH4201_2	12	46321441	T	G	SCAF11	17322	2203	3.81643E-05	0.00011398	1.2221E-08	8.6509E-05	0.0002149	0.132397333
30SCAF11_12_46321441_T_G.PH4201_3	12	46321441	T	G	SCAF11	89848	13541	1.49858E-05	4.9336E-05	2.4084E-09	8.6509E-05	0.0002149	0.132397333
40LAMA3_18_21453038_C_T.PH4201_1	18	21453038	C	T	LAMA3	120826	18663	1.83113E-05	6.70941E-05	4.4557E-09	8.97691E-05	0.000222999	0.143225
41LAMA3_18_21453038_C_T.PH4201_2	18	21453038	C	T	LAMA3	80326	10141	3.60021E-05	0.000115684	1.3223E-08	8.97691E-05	0.000222999	0.143225
42LAMA3_18_21453038_C_T.PH4201_3	18	21453038	C	T	LAMA3	96768	14415	3.92469E-05	8.11081E-05	6.4963E-09	8.97691E-05	0.000222999	0.143225
43FLNA_X_153587777_G_C.PH4201_1	X	153587777	G	C	FLNA	95577	17749	1.93432E-05	5.30161E-05	2.782E-09	6.76114E-05	0.000167956	0.19288
44FLNA_X_153587777_G_C.PH4201_2	X	153587777	G	C	FLNA	103754	20188	4.62967E-05	9.9895E-05	9.8432E-09	6.76114E-05	0.000167956	0.19288
45FLNA_X_153587777_G_C.PH4201_3	X	153587777	G	C	FLNA	40366	8007	1.29684E-05	3.31758E-05	1.0887E-09	6.76114E-05	0.000167956	0.19288
46ZNF223_19_44571260_C_A.PH4201_1	19	44571260	C	A	ZNF223	55900	4	1.48408E-05	2.97638E-05	8.7625E-10	6.19546E-05	0.000153904	0.0605694
47ZNF223_19_44571260_C_A.PH4201_2	19	44571260	C	A	ZNF223	47820	5	1.77988E-05	6.33744E-05	3.9655E-09	6.19546E-05	0.000153904	0.0605694
48ZNF223_19_44571260_C_A.PH4201_3	19	44571260	C	A	ZNF223	46731	11	3.29939E-05	8.22195E-05	6.6734E-09	6.19546E-05	0.000153904	0.0605694
52SCAF11_12_46321441_T_G.PH4201_1	12	46321441	T	G	SCAF11	70921	4630	1.46677E-05	3.74441E-05	1.385E-09	6.34898E-05	0.000157117	0.0605694
53SCAF11_12_46321441_T_G.PH4201_2	12	46321441	T	G	SCAF11	74343	4407	1.16077E-05	9.06887E-05	8.1176E-09	6.34898E-05	0.000157117	0.0605694
54SCAF11_12_46321441_T_G.PH4201_3	12	46321441	T	G	SCAF11	132680	7582	1.65308E-05	5.11932E-05	2.5903E-09	6.34898E-05	0.000157117	0.0605694
64LAMA3_18_21453038_C_T.PH4201_1	18	21453038	C	T	LAMA3	136420	9410	1.60786E-05	4.93645E-05	2.412E-09	7.02577E-05	0.00017453	0.0648923
65LAMA3_18_21453038_C_T.PH4201_2	18	21453038	C	T	LAMA3	95136	5138	2.3178E-05	7.34881E-05	5.3354E-09	7.02577E-05	0.00017453	0.0648923
66LAMA3_18_21453038_C_T.PH4201_3	18	21453038	C	T	LAMA3	103917	7450	3.96786E-05	8.45598E-05	7.061E-09	7.02577E-05	0.00017453	0.0648923
67FLNA_X_153587777_G_C.PH4201_1	X	153587777	G	C	FLNA	102340	8898	2.23604E-05	5.06317E-05	2.5374E-09	7.2709E-05	0.000180619	0.087697467
68FLNA_X_153587777_G_C.PH4201_2	X	153587777	G	C	FLNA	89619	8593	5.4927E-05	0.000109454	1.1807E-08	7.2709E-05	0.000180619	0.087697467
69FLNA_X_153587777_G_C.PH4201_3	X	153587777	G	C	FLNA	64276	5159	1.58968E-05	3.91355E-05	1.5158E-09	7.2709E-05	0.000180619	0.087697467
70ZNF223_19_44571260_C_A.PH4201_1	19	44571260	C	A	ZNF223	44181	2	1.44693E-05	2.98906E-05	8.8363E-10	5.28463E-05	0.000131278	3.4157E-05
71ZNF223_19_44571260_C_A.PH4201_2	19	44571260	C	A	ZNF223	52445	3	1.4025E-05	4.72239E-05	2.2019E-09	5.28463E-05	0.000131278	3.4157E-05
72ZNF223_19_44571260_C_A.PH4201_3	19	44571260	C	A	ZNF223	45498	0	2.60957E-05	7.32218E-05	5.2927E-09	5.28463E-05	0.000131278	3.4157E-05
73FLNA_X_153579448_A_G.PH4201_1	X	153579448	A	G	FLNA	107654	6311	1.06204E-05	8.28234E-05	7.8877E-10	3.45601E-05	8.58522E-05	0.045955333
74FLNA_X_153579448_A_G.PH4201_2	X	153579448	A	G	FLNA	103932	4578	1.35519E-05	8.30512E-05	1.4316E-09	3.45601E-05	8.58522E-05	0.045955333
75FLNA_X_153579448_A_G.PH4201_3	X	153579448	A	G	FLNA	57423	2021	1.2331E-05	3.71279E-05	1.3628E-09	3.45601E-05	8.58522E-05	0.045955333
82SLX4_16_3639306_G_A.PH4201_1	16	3639306	G	A	SLX4	114402	3794	4.05151E-05	6.0057E-05	3.5689E-09	9.85036E-05	0.000244696	0.031055833
83SLX4_16_3639306_G_A.PH4201_2	16	3639306	G	A	SLX4	122229	3343	2.16981E-05	6.56154E-05	4.2509E-09	9.85036E-05	0.000244696	0.031055833
84SLX4_16_3639306_G_A.PH4201_3	16	3639306	G	A	SLX4	107799	3520	4.92961E-05	0.00014669	2.1289E-08	9.85036E-05	0.000244696	0.031055833
88LAMA3_18_21453038_C_T.PH4201_1	18	21453038	C	T	LAMA3	141739	3518	1.6245E-05	4.75211E-05	2.2352E-09	6.60986E-05	0.000164198	0.031337033
89LAMA3_18_21453038_C_T.PH4201_2	18	21453038	C	T	LAMA3	123130	4064	1.96342E-05	5.26053E-05	2.7372E-09	6.60986E-05	0.000164198	0.031337033
90LAMA3_18_21453038_C_T.PH4201_3	18	21453038	C	T	LAMA3	96504	3492	4.10938E-05	9.07612E-05	8.1346E-09	6.60986E-05	0.000164198	0.031337033
91FLNA_X_153587777_G_C.PH4201_1	X	153587777	G	C	FLNA	120137	5731	1.54135E-05	3.71281E-05	1.3644E-09	5.67267E-05	0.000140917	0.0437276
92FLNA_X_153587777_G_C.PH4201_2	X	153587777	G	C	FLNA	144879	6360	3.3529E-05	8.33444E-05	6.8511E-09	5.67267E-05	0.000140917	0.0437276
93FLNA_X_153587777_G_C.PH4201_3	X	153587777	G	C	FLNA	78221	3096	2.01674E-05	3.81206E-05	1.4382E-09	5.67267E-05	0.000140917	0.0437276
94ZNF223_19_44571260_C_A.PH4201_1	19	44571260	C	A	ZNF223	42947	2	1.60091E-05	4.1851E-05	1.7333E-09	5.59136E-05	0.000138897	3.31067E-05
95ZNF223_19_44571260_C_A.PH4201_2	19	44571260	C	A	ZNF223	65985	0	1.81489E-05	5.1892E-05	2.6559E-09	5.59136E-05	0.000138897	3.31067E-05
96ZNF223_19_44571260_C_A.PH4201_3	19	44571260	C	A	ZNF223	56871	3	3.58294E-05	7.11016E-05	4.9898E-09	5.59136E-05	0.000138897	3.31067E-05
97FLNA_X_153579448_A_G.PH4201_1	X	153579448	A	G	FLNA	112175	4374	7.44648E-06	2.09631E-05	4.3492E-10	2.97837E-05	7.39868E-05	0.030000867
98FLNA_X_153579448_A_G.PH4201_2	X	153579448	A	G	FLNA	112537	3127	1.3643E-05	3.99966E-05	1.5827E-09	2.97837E-05	7.39868E-05	0.030000867
99FLNA_X_153579448_A_G.PH4201_3	X	153579448	A	G	FLNA	74364	1727	1.24087E-05	2.55048E-05	6.4357E-10	2.97837E-05	7.39868E-05	0.030000867
100SCAF11_12_46321441_T_G.PH4201_1	12	46321441	T	G	SCAF11	105758	1337	6.32158E-06	1.39944E-05	1.9364E-10	7.36666E-05	0.000181507	0.012333833

TABLE 5-continued

Observed alternate allele fractions												
PrimerID	Chr	AlleleStart	Ref	Alt Gene	IT Read Depth	Alt Allele Depth	Background AAF (within 50 nts)	Stdev Background AAF	Variance Background AAF	Stdev of Average Background	Confidence interval of Background	Average AAF
101SCAF11_12_46321441_T_G_P44201_2	12	46321441	T	G SCAF11	17613	203	3.59949E-05	0.00011086	1.213E-08	7.30666E-05	0.000181507	0.012333833
102SCAF11_12_46321441_T_G_P44201_3	12	46321441	T	G SCAF11	139164	1786	2.03063E-05	6.11195E-05	3.6922E-09	7.30666E-05	0.000181507	0.012333833
112SCAF11_18_21453038_C_T_P44201_1	18	21453038	C	T LAMA3	144251	2632	1.19654E-05	4.35623E-05	1.8783E-09	0.000123515	0.000306828	0.011833627
113LAMA3_18_21453038_C_T_P44201_2	18	21453038	C	T LAMA3	13958	18	0.000131449	0.000195468	3.7537E-08	0.000123515	0.000306828	0.011833627
114LAMA3_18_21453038_C_T_P44201_3	18	21453038	C	T LAMA3	126650	2022	3.34354E-05	8.02018E-05	6.3519E-09	0.000123515	0.000306828	0.011833627
115FLNA_X_153587777_G_C_P44201_1	X	153587777	G	C FLNA	117886	2782	1.16224E-05	2.79385E-05	7.726E-10	8.47043E-05	0.000210417	0.020187733
117FLNA_X_153587777_G_C_P44201_3	X	153587777	G	C FLNA	82453	1404	1.44532E-05	3.62043E-05	1.2972E-09	8.47043E-05	0.000210417	0.020187733
118ZNF223_19_44571260_C_A_P44201_1	19	44571260	C	A ZNF223	46424	0	1.75212E-05	3.98215E-05	1.5687E-09	6.15054E-05	0.000152788	2.88507E-05
119ZNF223_19_44571260_C_A_P44201_2	19	44571260	C	A ZNF223	48037	0	1.41639E-05	5.2454E-05	2.7166E-09	6.15054E-05	0.000152788	2.88507E-05
120ZNF223_19_44571260_C_A_P44201_3	19	44571260	C	A ZNF223	46215	4	3.01323E-05	8.45884E-05	7.0635E-09	6.15054E-05	0.000152788	2.88507E-05
121FLNA_X_153579448_A_G_P44201_1	X	153579448	A	G FLNA	104910	1543	1.23152E-05	3.51172E-05	1.2241E-09	4.31322E-05	0.000107146	0.011726387
122FLNA_X_153579448_A_G_P44201_2	X	153579448	A	G FLNA	87067	1273	2.03354E-05	6.03457E-05	3.6011E-09	4.31322E-05	0.000107146	0.011726387
123FLNA_X_153579448_A_G_P44201_3	X	153579448	A	G FLNA	60679	355	1.4749E-05	2.76533E-05	7.5592E-10	4.31322E-05	0.000107146	0.011726387
124SCAF11_12_46321441_T_G_P44201_1	12	46321441	T	G SCAF11	128949	708	6.19744E-05	0.000144225	2.0567E-08	9.41463E-05	0.000233872	0.006252057
125SCAF11_12_46321441_T_G_P44201_2	12	46321441	T	G SCAF11	86647	600	2.55061E-05	6.46009E-05	4.1191E-09	9.41463E-05	0.000233872	0.006252057
126SCAF11_12_46321441_T_G_P44201_3	12	46321441	T	G SCAF11	138622	879	1.41111E-05	4.38961E-05	1.9045E-09	9.41463E-05	0.000233872	0.006252057
130SLX4_16_3639306_G_A_P44201_1	16	3639306	G	A SLX4	90970	706	3.60242E-05	6.68836E-05	4.4258E-09	7.93998E-05	0.00019724	0.00776187
131SLX4_16_3639306_G_A_P44201_2	16	3639306	G	A SLX4	107933	865	1.41706E-05	5.09717E-05	2.5652E-09	7.93998E-05	0.00019724	0.00776187
132SLX4_16_3639306_G_A_P44201_3	16	3639306	G	A SLX4	124491	935	4.18812E-05	0.000109779	1.1922E-08	7.93998E-05	0.00019724	0.00776187
136LAMA3_18_21453038_C_T_P44201_1	18	21453038	C	T LAMA3	145645	1114	1.19383E-05	4.08667E-05	1.653E-09	7.02583E-05	0.000174531	0.007249317
137LAMA3_18_21453038_C_T_P44201_2	18	21453038	C	T LAMA3	94621	562	3.00284E-05	6.81186E-05	4.5849E-09	7.02583E-05	0.000174531	0.007249317
138LAMA3_18_21453038_C_T_P44201_3	18	21453038	C	T LAMA3	56742	463	3.56723E-05	9.31623E-05	8.5707E-09	7.02583E-05	0.000174531	0.007249317
139FLNA_X_153587777_G_C_P44201_1	X	153587777	G	C FLNA	78250	822	1.89315E-05	3.97105E-05	1.5608E-09	4.86093E-05	0.000120752	0.01034064
140FLNA_X_153587777_G_C_P44201_2	X	153587777	G	C FLNA	135680	1295	2.90505E-05	6.76937E-05	4.5197E-09	4.86093E-05	0.000120752	0.01034064
141FLNA_X_153587777_G_C_P44201_3	X	153587777	G	C FLNA	90316	991	1.53494E-05	3.19139E-05	1.0081E-09	4.86093E-05	0.000120752	0.01034064
142ZNF223_19_44571260_C_A_P44201_1	19	44571260	C	A ZNF223	45173	140	3.13548E-05	4.1112E-05	1.6728E-09	5.27123E-05	0.000130945	0.001166528
143ZNF223_19_44571260_C_A_P44201_2	19	44571260	C	A ZNF223	66182	2	1.31595E-05	4.74401E-05	2.2221E-09	5.27123E-05	0.000130945	0.001166528
144ZNF223_19_44571260_C_A_P44201_3	19	44571260	C	A ZNF223	32418	12	2.60104E-05	6.71356E-05	4.4409E-09	5.27123E-05	0.000130945	0.001166528
182NA_5_73717969_G_A_C-17_3	5	73717969	G	A NA	133924	9487	1.61668E-05	5.32435E-05	2.8056E-09	6.87898E-05	0.000170883	0.0728698
183NA_5_73717969_G_A_C-18_3	5	73717969	G	A NA	148542	15441	1.06255E-05	4.24466E-05	1.7833E-09	5.751E-05	0.000142863	0.105353
184NA_5_73717969_G_A_C-9_3	5	73717969	G	A NA	149125	16289	1.72387E-05	4.24724E-05	1.7855E-09	9.5883E-05	0.000238192	0.1056935
185NA_5_73717969_G_A_C-11_3	5	73717969	G	A NA	149150	16863	1.20897E-05	4.73091E-05	2.2153E-09	5.2303E-05	0.000129928	0.1146715
187NA_5_73717969_G_A_C-45_3	5	73717969	G	A NA	148515	18657	1.03337E-05	4.0947E-05	1.6596E-09	5.2003E-05	0.000129184	0.122451
189NA_5_73717969_G_A_C-17_1	5	73717969	G	A NA	127236	8870	3.7739E-05	0.000106012	1.1096E-08	6.87898E-05	0.000170883	0.0728698
190NA_5_73717969_G_A_C-18_1	5	73717969	G	A NA	128246	13691	3.03114E-05	6.99529E-05	4.8315E-09	5.751E-05	0.000142863	0.105353
191NA_5_73717969_G_A_C-9_1	5	73717969	G	A NA	126424	12915	4.73437E-05	0.000129674	1.6602E-08	9.5883E-05	0.000238192	0.1056935
192NA_5_73717969_G_A_C-11_1	5	73717969	G	A NA	129169	15020	2.79676E-05	5.7425E-05	3.2559E-09	5.2303E-05	0.000129928	0.1146715
194NA_5_73717969_G_A_C-45_1	5	73717969	G	A NA	127861	15251	2.9291E-05	6.16219E-05	3.7492E-09	5.2003E-05	0.000129184	0.122451
196NA_5_73717969_G_A_C-17_2	5	73717969	G	A NA	146571	11441	6.38217E-06	1.72436E-05	2.9425E-10	6.87898E-05	0.000170883	0.0728698
199NA_11_49854989_C_T_A-9_3	11	49854989	C	T NA	32775	3	1.75167E-05	4.59566E-05	2.0895E-09	3.59978E-05	8.94235E-05	4.57666E-05
200NA_11_49854989_C_T_A-9_1	11	49854989	C	T NA	141978	0	1.14214E-05	2.25275E-05	5.0215E-10	3.59978E-05	8.94235E-05	4.57666E-05
204NA_11_170130646_T_G_C-9_3	11	170130646	T	G NA	147507	0	1.36308E-05	3.00914E-05	8.896E-10	2.98262E-05	7.40925E-05	0
1FLNA_X_153579448_A_G_P44201_1	X	153579448	A	G FLNA	24264	287	6.34352E-05	1.3124E-09	3.04522E-05	7.56476E-05	0.007530833	0.007530833
2FLNA_X_153579448_A_G_P44201_2	X	153579448	A	G FLNA	69733	543	5.2092E-06	1.02847E-05	1.0461E-10	3.04522E-05	7.56476E-05	0.007530833
3FLNA_X_153579448_A_G_P44201_3	X	153579448	A	G FLNA	65828	196	1.1365E-05	3.71572E-05	1.365E-09	3.04522E-05	7.56476E-05	0.007530833
4SCAF11_12_46321441_T_G_P44201_1	12	46321441	T	G SCAF11	81926	249	3.69937E-06	9.58745E-06	9.0886E-11	8.90779E-05	0.000221282	0.003211453
5SCAF11_12_46321441_T_G_P44201_2	12	46321441	T	G SCAF11	15296	51	3.52443E-05	9.67775E-05	9.2426E-09	8.90779E-05	0.000221282	0.003211453

TABLE 5-continued

Observed alternate allele fractions														
PrimerID	Chr	AlleleStart	Ref	Alt	Gene	IT Read Depth	Alt Allele Depth	Background AAF (within 50 nts)	Stdev AAF	Variance Background AAF	Stdev of Average Background	Confidence interval of		Average AAF
												Background	Background	
68SCAF11_12_46321441_T_G_P44201_3	12	46321441	T	G	SCAF11	133402	435	2.51523E-05	0.000120985	1.4471E-08	8.90779E-05	0.000221282	0.003211453	
108SLX4_16_3639306_G_A_P44201_1	16	3639306	G	A	SLX4	58364	322	2.67548E-05	7.55218E-05	5.6429E-09	6.73416E-05	0.000167286	0.00393089	
118SLX4_16_3639306_G_A_P44201_2	16	3639306	G	A	SLX4	44842	110	2.6859E-05	6.34548E-05	3.9749E-09	6.73416E-05	0.000167286	0.00393089	
12SLX4_16_3639306_G_A_P44201_3	16	3639306	G	A	SLX4	59385	227	3.10302E-05	6.34771E-05	3.9869E-09	6.73416E-05	0.000167286	0.00393089	
16LAMA3_18_21453038_C_T_P44201_1	18	21453038	C	T	LAMA3	95634	485	2.04848E-05	7.29925E-05	5.2735E-09	6.74429E-05	0.000167537	0.004350827	
17LAMA3_18_21453038_C_T_P44201_2	18	21453038	C	T	LAMA3	37744	145	2.20973E-05	5.66442E-05	3.1614E-09	6.74429E-05	0.000167537	0.004350827	
18LAMA3_18_21453038_C_T_P44201_3	18	21453038	C	T	LAMA3	75615	313	2.89559E-05	7.26407E-05	5.2107E-09	6.74429E-05	0.000167537	0.004350827	
19FLNA_X_153587777_G_C_P44201_1	X	153587777	G	C	FLNA	91142	673	9.17142E-06	2.32548E-05	5.3527E-10	9.82176E-05	0.000243986	0.006361423	
20FLNA_X_153587777_G_C_P44201_2	X	153587777	G	C	FLNA	99527	756	5.65447E-05	0.000164486	2.668E-08	9.82176E-05	0.000243986	0.006361423	
21FLNA_X_153587777_G_C_P44201_3	X	153587777	G	C	FLNA	79186	325	1.4775E-05	4.17448E-05	1.7248E-09	9.82176E-05	0.000243986	0.006361423	
22ZNF223_19_44571260_C_A_P44201_1	19	44571260	C	A	ZNF223	46605	1	8.77945E-06	2.22494E-05	4.8988E-10	5.00911E-05	0.000124433	1.26338E-05	
23ZNF223_19_44571260_C_A_P44201_2	19	44571260	C	A	ZNF223	70797	0	1.24727E-05	4.14129E-05	1.693E-09	5.00911E-05	0.000124433	1.26338E-05	
24ZNF223_19_44571260_C_A_P44201_3	19	44571260	C	A	ZNF223	60811	1	2.99432E-05	7.36047E-05	5.3444E-09	5.00911E-05	0.000124433	1.26338E-05	
25FLNA_X_153579448_A_G_P44201_1	X	153579448	A	G	FLNA	129429	464	8.16419E-06	2.69819E-05	7.2052E-10	2.30552E-05	5.72723E-05	0.00236643	
26FLNA_X_153579448_A_G_P44201_2	X	153579448	A	G	FLNA	95869	199	1.05279E-05	1.80926E-05	3.2386E-10	2.30552E-05	5.72723E-05	0.00236643	
27FLNA_X_153579448_A_G_P44201_3	X	153579448	A	G	FLNA	59087	85	1.40795E-05	2.35888E-05	5.5025E-10	2.30552E-05	5.72723E-05	0.00236643	
28SCAF11_12_46321441_T_G_P44201_1	12	46321441	T	G	SCAF11	83268	112	4.5364E-06	1.50852E-05	2.2501E-10	4.26917E-05	0.000106052	0.001588223	
30SCAF11_12_46321441_T_G_P44201_3	12	46321441	T	G	SCAF11	117362	122	1.16431E-05	3.55168E-05	1.2473E-09	4.26917E-05	0.000106052	0.001588223	
34SLX4_16_3639306_G_A_P44201_1	16	3639306	G	A	SLX4	88998	170	2.1575E-05	6.08845E-05	3.6657E-09	8.8215E-05	0.000219138	0.001983499	
35SLX4_16_3639306_G_A_P44201_2	16	3639306	G	A	SLX4	67120	56	4.66127E-05	0.000122221	1.4746E-08	8.8215E-05	0.000219138	0.001983499	
36SLX4_16_3639306_G_A_P44201_3	16	3639306	G	A	SLX4	61759	198	3.13336E-05	7.0612E-05	4.9336E-09	8.8215E-05	0.000219138	0.001983499	
43FLNA_X_153587777_G_C_P44201_1	X	153587777	G	C	FLNA	107980	337	1.79554E-05	3.73283E-05	1.3792E-09	4.81874E-05	0.000119704	0.00332788	
44FLNA_X_153587777_G_C_P44201_2	X	153587777	G	C	FLNA	132263	411	2.41059E-05	6.49675E-05	4.163E-09	4.81874E-05	0.000119704	0.00332788	
45FLNA_X_153587777_G_C_P44201_3	X	153587777	G	C	FLNA	68704	258	1.47833E-05	3.79376E-05	1.424E-09	4.81874E-05	0.000119704	0.00332788	
46ZNF223_19_44571260_C_A_P44201_1	19	44571260	C	A	ZNF223	57574	0	1.03416E-05	2.74209E-05	7.4407E-10	5.19426E-05	0.000129032	6.7653E-06	
47ZNF223_19_44571260_C_A_P44201_2	19	44571260	C	A	ZNF223	69787	0	1.12838E-05	4.2106E-05	1.7502E-09	5.19426E-05	0.000129032	6.7653E-06	
48ZNF223_19_44571260_C_A_P44201_3	19	44571260	C	A	ZNF223	49271	1	3.3871E-05	7.53571E-05	5.5998E-09	5.19426E-05	0.000129032	6.7653E-06	
49FLNA_X_153579448_A_G_P44201_1	X	153579448	A	G	FLNA	87941	197	1.1037E-05	2.80666E-05	7.7953E-10	2.67942E-05	0.000149334	0.001935583	
50FLNA_X_153579448_A_G_P44201_2	X	153579448	A	G	FLNA	85699	173	1.1882E-05	2.72096E-05	7.3257E-10	2.67942E-05	0.000149334	0.001935583	
51FLNA_X_153579448_A_G_P44201_3	X	153579448	A	G	FLNA	52298	81	1.08666E-05	2.54719E-05	6.4169E-10	2.67942E-05	0.000149334	0.001935583	
52SCAF11_12_46321441_T_G_P44201_1	12	46321441	T	G	SCAF11	54590	56	1.5484E-05	3.75274E-05	1.3907E-09	6.01151E-05	0.000149334	0.000855501	
53SCAF11_12_46321441_T_G_P44201_2	12	46321441	T	G	SCAF11	35287	30	3.05225E-05	8.63868E-05	7.3638E-09	6.01151E-05	0.000149334	0.000855501	
54SCAF11_12_46321441_T_G_P44201_3	12	46321441	T	G	SCAF11	130340	90	1.40915E-05	4.59206E-05	2.085E-09	6.01151E-05	0.000149334	0.000855501	
67FLNA_X_153587777_G_C_P44201_1	X	153587777	G	C	FLNA	119586	240	1.18964E-05	3.5427E-05	1.2423E-09	6.63438E-05	0.000164807	0.00198682	
68FLNA_X_153587777_G_C_P44201_2	X	153587777	G	C	FLNA	95355	226	4.46509E-05	0.000105779	1.1032E-08	6.63438E-05	0.000164807	0.00198682	
69FLNA_X_153587777_G_C_P44201_3	X	153587777	G	C	FLNA	56838	90	1.16134E-05	3.06666E-05	9.3065E-10	6.63438E-05	0.000164807	0.00198682	
70ZNF223_19_44571260_C_A_P44201_1	19	44571260	C	A	ZNF223	59122	0	1.20736E-05	2.8092E-05	7.8094E-10	4.38919E-05	0.000109033	0	
71ZNF223_19_44571260_C_A_P44201_2	19	44571260	C	A	ZNF223	64656	0	1.0793E-05	3.9434E-05	1.5354E-09	4.38919E-05	0.000109033	0	
72ZNF223_19_44571260_C_A_P44201_3	19	44571260	C	A	ZNF223	58814	0	2.39246E-05	5.92561E-05	3.4632E-09	4.38919E-05	0.000109033	0	
73FLNA_X_153579448_A_G_P44201_1	X	153579448	A	G	FLNA	86773	167	9.98466E-06	2.56786E-05	6.5252E-10	3.30698E-05	8.21499E-05	0.001185302	
74FLNA_X_153579448_A_G_P44201_2	X	153579448	A	G	FLNA	72101	70	1.20489E-05	3.55212E-05	1.2474E-09	3.30698E-05	8.21499E-05	0.001185302	
75FLNA_X_153579448_A_G_P44201_3	X	153579448	A	G	FLNA	42393	28	1.63782E-05	3.73709E-05	1.3809E-09	3.30698E-05	8.21499E-05	0.001185302	
88LAMA3_18_21453038_C_T_P44201_1	18	21453038	C	T	LAMA3	82365	74	1.89137E-05	5.99713E-05	3.5599E-09	6.7586E-05	0.000167893	0.00088785	
89LAMA3_18_21453038_C_T_P44201_2	18	21453038	C	T	LAMA3	73739	73	2.10443E-05	5.8422E-05	3.3748E-09	6.7586E-05	0.000167893	0.00088785	
90LAMA3_18_21453038_C_T_P44201_3	18	21453038	C	T	LAMA3	98048	76	3.09949E-05	8.27928E-05	6.769E-09	6.7586E-05	0.000167893	0.00088785	
91FLNA_X_153587777_G_C_P44201_1	X	153587777	G	C	FLNA	96939	4	1.70301E-05	4.40317E-05	1.919E-09	6.08077E-05	0.000151055	0.000592997	
92FLNA_X_153587777_G_C_P44201_2	X	153587777	G	C	FLNA	108834	94	3.7497E-05	9.10425E-05	8.1752E-09	6.08077E-05	0.000151055	0.000592997	

TABLE 5-continued

Observed alternate allele fractions													
PrimerID	Chr	AlleleStart	Ref	Alt Gene	IT Read Depth	Alt Allele Depth	Background AAF (within 50 nts)	Stdev Background AAF	Variance Background AAF	Stdev of Average		Confidence interval of Background	Average AAF
										Background	Background		
93FLNA_X_153587777_G_C_P44201_3	X	153587777	G	C FLNA	69792	61	1.29442E-05	3.17658E-05	9.9855E-10	6.08077E-05	0.000151055	0.000592997	0
94ZNF223_19_44571260_C_A_P44201_1	19	44571260	C	A ZNF223	59496	0	1.31575E-05	2.61998E-05	6.7928E-10	5.65725E-05	0.000140534	0	0
95ZNF223_19_44571260_C_A_P44201_2	19	44571260	C	A ZNF223	60176	0	1.09985E-05	4.22833E-05	1.7644E-09	5.65725E-05	0.000140534	0	0
96ZNF223_19_44571260_C_A_P44201_3	19	44571260	C	A ZNF223	53680	0	3.38316E-05	8.51887E-05	7.1577E-09	5.65725E-05	0.000140534	0	0
97FLNA_X_153579448_A_G_P44201_1	X	153579448	A	G FLNA	97371	77	9.51036E-06	2.66569E-05	7.0327E-10	3.44963E-05	8.56935E-05	0.000477369	0
98FLNA_X_153579448_A_G_P44201_2	X	153579448	A	G FLNA	63812	24	1.02225E-05	2.64331E-05	6.9127E-10	3.44963E-05	8.56935E-05	0.000477369	0
99FLNA_X_153579448_A_G_P44201_3	X	153579448	A	G FLNA	26394	7	2.43236E-05	4.69184E-05	2.1754E-09	3.44963E-05	8.56935E-05	0.000477369	0
102SCAF11_12_46321441_T_G_P44201_3	12	46321441	T	G SCAF11	28326	0	1.8436E-05	5.25987E-05	2.7355E-09	6.88974E-05	0.000166182	0.000250532	0
106SLX4_16_3639306_G_A_P44201_1	16	3639306	G	A SLX4	69218	6	1.20893E-05	3.92501E-05	1.5236E-09	4.88915E-05	0.000121453	0.000102273	0
107SLX4_16_3639306_G_A_P44201_2	16	3639306	G	A SLX4	38925	0	1.73624E-05	5.98804E-05	3.5403E-09	4.88915E-05	0.000121453	0.000102273	0
108SLX4_16_3639306_G_A_P44201_3	16	3639306	G	A SLX4	54512	12	1.79625E-05	4.61508E-05	2.1072E-09	4.88915E-05	0.000121453	0.000102273	0
115FLNA_X_153587777_G_C_P44201_1	X	153587777	G	C FLNA	80552	1	6.90677E-06	2.10375E-05	4.3806E-10	6.01185E-05	0.000149343	0.000400965	0
117FLNA_X_153587777_G_C_P44201_3	X	153587777	G	C FLNA	58499	0	9.26507E-06	1.90282E-05	3.583E-10	6.01185E-05	0.000149343	0.000400965	0
119ZNF223_19_44571260_C_A_P44201_2	19	44571260	C	A ZNF223	41307	0	1.1469E-05	4.05857E-05	1.6261E-09	7.439E-05	0.000184795	0	0
121FLNA_X_153579448_A_G_P44201_1	X	153579448	A	G FLNA	83614	23	7.52473E-06	1.95783E-05	3.7928E-10	1.73495E-05	4.30987E-05	0.000191769	0
122FLNA_X_153579448_A_G_P44201_2	X	153579448	A	G FLNA	58905	9	6.48471E-06	1.69247E-05	2.8319E-10	1.73495E-05	4.30987E-05	0.000191769	0
123FLNA_X_153579448_A_G_P44201_3	X	153579448	A	G FLNA	54258	8	5.38074E-06	1.55977E-05	2.4055E-10	1.73495E-05	4.30987E-05	0.000191769	0
124SCAF11_12_46321441_T_G_P44201_1	12	46321441	T	G SCAF11	10183	0	1.14584E-05	3.89589E-05	1.5007E-09	5.42166E-05	0.000134681	4.21327E-05	0
125SCAF11_12_46321441_T_G_P44201_2	12	46321441	T	G SCAF11	15823	2	2.08904E-05	7.62898E-05	5.7459E-09	5.42166E-05	0.000134681	4.21327E-05	0
126SCAF11_12_46321441_T_G_P44201_3	12	46321441	T	G SCAF11	43369	0	1.05994E-05	3.98686E-05	1.5716E-09	5.42166E-05	0.000134681	4.21327E-05	0
130SLX4_16_3639306_G_A_P44201_1	16	3639306	G	A SLX4	42958	8	1.40447E-05	3.79384E-05	1.4238E-09	3.76729E-05	9.35848E-05	0.000135465	0
131SLX4_16_3639306_G_A_P44201_2	16	3639306	G	A SLX4	22424	2	1.55039E-05	3.99029E-05	1.5721E-09	3.76729E-05	9.35848E-05	0.000135465	0
132SLX4_16_3639306_G_A_P44201_3	16	3639306	G	A SLX4	53444	7	1.75157E-05	3.57126E-05	1.2618E-09	3.76729E-05	9.35848E-05	0.000135465	0
136LAMA3_18_21453038_C_T_P44201_1	18	21453038	C	T LAMA3	121601	80	9.09633E-06	2.94624E-05	8.5917E-10	5.14128E-05	0.000127716	0.000420001	0
137LAMA3_18_21453038_C_T_P44201_2	18	21453038	C	T LAMA3	53584	12	2.4705E-05	6.33785E-05	3.9684E-09	5.14128E-05	0.000127716	0.000420001	0
138LAMA3_18_21453038_C_T_P44201_3	18	21453038	C	T LAMA3	47598	18	2.24107E-05	5.60533E-05	3.1022E-09	5.14128E-05	0.000127716	0.000420001	0
139FLNA_X_153587777_G_C_P44201_1	X	153587777	G	C FLNA	71859	2	9.87313E-06	2.39318E-05	5.6689E-10	5.99534E-05	0.000148933	0.000210718	0
140FLNA_X_153587777_G_C_P44201_2	X	153587777	G	C FLNA	54664	31	3.24141E-05	9.86693E-05	9.6023E-09	5.99534E-05	0.000148933	0.000210718	0
141FLNA_X_153587777_G_C_P44201_3	X	153587777	G	C FLNA	53732	2	9.80512E-06	2.49083E-05	6.1409E-10	5.99534E-05	0.000148933	0.000210718	0
142ZNF223_19_44571260_C_A_P44201_1	19	44571260	C	A ZNF223	31245	1	1.11919E-05	2.51831E-05	6.2765E-10	0.000189244	0.000470108	0.000419128	0
143ZNF223_19_44571260_C_A_P44201_2	19	44571260	C	A ZNF223	37757	0	1.22571E-05	5.98359E-05	3.3544E-09	0.000189244	0.000470108	0.000419128	0
144ZNF223_19_44571260_C_A_P44201_3	19	44571260	C	A ZNF223	11425	14	0.000119374	0.000323909	1.0328E-07	0.000189244	0.000470108	0.000419128	0
182NA_5_73717969_G_A_C-putamen_3	5	73717969	G	A NA	146868	16657	1.04913E-05	3.63854E-05	1.3103E-09	7.32249E-05	0.000181901	0.1121115	0
183NA_5_73717969_G_A_C-37_3	5	73717969	G	A NA	148025	16387	9.50264E-06	3.4586E-05	1.184E-09	4.9362E-05	0.000122622	0.1104095	0
184NA_5_73717969_G_A_C-7_3	5	73717969	G	A NA	148518	14328	1.04058E-05	3.78078E-05	1.4148E-09	6.39402E-05	0.000158836	0.0963384	0
185NA_5_73717969_G_A_C-19_3	5	73717969	G	A NA	148641	19027	1.21165E-05	4.49393E-05	1.9989E-09	5.8551E-05	0.000145449	0.1281165	0
186NA_5_73717969_G_A_C-pois_3	5	73717969	G	A NA	146034	15926	9.58084E-06	3.91421E-05	1.5165E-09	0.000106587	0.000264777	0.1114735	0
187NA_5_73717969_G_A_C-adrenal_3	5	73717969	G	A NA	148408	18181	9.66387E-06	4.03626E-05	1.6125E-09	4.5798E-05	0.000113768	0.1167425	0
188NA_5_73717969_G_A_C-pancreas_3	5	73717969	G	A NA	139170	14360	1.77384E-05	5.30139E-05	2.7812E-09	0.000144591	0.000359185	0.09964535	0
189NA_5_73717969_G_A_C-putamen_1	5	73717969	G	A NA	131453	14566	4.26031E-05	9.76354E-05	9.4135E-09	7.32249E-05	0.000181901	0.1121115	0
190NA_5_73717969_G_A_C-37_1	5	73717969	G	A NA	135104	14877	2.94044E-05	6.11223E-05	3.6892E-09	4.9362E-05	0.000122622	0.1104095	0
191NA_5_73717969_G_A_C-7_1	5	73717969	G	A NA	132828	12726	3.32702E-05	8.27493E-05	6.7619E-09	6.39402E-05	0.000158836	0.0963384	0
192NA_5_73717969_G_A_C-19_1	5	73717969	G	A NA	134589	17258	2.72636E-05	7.01355E-05	4.8575E-09	5.8551E-05	0.000145449	0.1281165	0
193NA_5_73717969_G_A_C-pois_1	5	73717969	G	A NA	130810	14898	4.3965E-05	0.000146539	2.1205E-08	0.000106587	0.000264777	0.1114735	0
194NA_5_73717969_G_A_C-adrenal_1	5	73717969	G	A NA	134856	14966	2.2429E-05	5.111379E-05	2.5824E-09	4.5798E-05	0.000113768	0.1167425	0
195NA_5_73717969_G_A_C-pancreas_1	5	73717969	G	A NA	130978	12588	5.25575E-05	0.000198812	3.9032E-08	0.000144591	0.000359185	0.09964535	0
196NA_5_73717969_G_A_C-cerebellum_2	5	73717969	G	A NA	145082	21160	8.00971E-06	2.02271E-05	4.0488E-10	2.01215E-05	4.99846E-05	0.145849	0

TABLE 5-continued

PrimerID	Chr	AlleleStart	Ref	Alt Gene	IT Read Depth	Alt Allele Depth	Background AAF (within 50 ms)	Stdev Background AAF	Variance Background AAF	Stdev of Average Background	Confidence interval of Background	Average AAF
200NA_3_177844577_G_A_C-17_1	3	177844577	G	A NA	84017	46	1.38304E-05	4.52855E-05	2.0272E-09	3.85259E-05	9.57036E-05	0.000361819
201NA_3_177844577_G_A_C-18_1	3	177844577	G	A NA	146570	343	6.03893E-06	1.68273E-05	2.8015E-10	1.33039E-05	3.30488E-05	0.00252737
202NA_3_177844577_G_A_C-19_1	3	177844577	G	A NA	140940	1026	6.91288E-06	2.74597E-05	7.4597E-10	2.89933E-05	7.20232E-05	0.007887165
203NA_3_177844577_G_A_C-11_1	3	177844577	G	A NA	133617	1015	1.17129E-05	3.84869E-05	1.4612E-09	2.82699E-05	7.02263E-05	0.01191867
204NA_3_177844577_G_A_C-47_1	3	177844577	G	A NA	146353	1231	8.90321E-06	2.56134E-05	6.4921E-10	5.15293E-05	0.000128006	0.010443985
205NA_3_177844577_G_A_C-45_1	3	177844577	G	A NA	136805	185	6.38948E-06	3.04637E-05	9.1761E-10	2.19555E-05	5.45392E-05	0.001199725
206NA_3_177844577_G_A_C-44_1	3	177844577	G	A NA	141268	395	1.73897E-06	5.76108E-06	3.2841E-11	1.20403E-05	2.99097E-05	0.002975415
207NA_3_177844577_G_A_C-17_3	3	177844577	G	A NA	45421	8	7.77866E-06	3.08413E-05	9.4128E-10	3.85259E-05	9.57036E-05	0.000361819
208NA_3_177844577_G_A_C-18_3	3	177844577	G	A NA	135197	367	2.28091E-06	8.63841E-06	7.3845E-11	1.33039E-05	3.30488E-05	0.00252737
209NA_3_177844577_G_A_C-9_3	3	177844577	G	A NA	60391	513	7.70605E-06	3.07423E-05	9.3524E-10	2.89933E-05	7.20232E-05	0.007887165
210NA_3_177844577_G_A_C-11_3	3	177844577	G	A NA	102580	1666	3.05559E-06	1.17725E-05	1.3715E-10	2.82699E-05	7.02263E-05	0.01191867
211NA_3_177844577_G_A_C-47_3	3	177844577	G	A NA	26449	330	3.77525E-05	6.86744E-05	4.6613E-09	5.15293E-05	0.000128006	0.010443985
212NA_3_177844577_G_A_C-45_3	3	177844577	G	A NA	121280	127	2.40403E-06	6.84998E-06	4.6433E-11	2.19555E-05	5.45392E-05	0.001199725
213NA_3_177844577_G_A_C-44_3	3	177844577	G	A NA	100801	318	3.94254E-06	1.61183E-05	2.5709E-10	1.20403E-05	2.99097E-05	0.002975415
215NA_3_177844577_G_A_C-8_3	3	177844577	G	A NA	17068	94	4.22828E-06	1.8048E-05	3.2234E-10	4.47728E-05	0.000112222	0.0051468
184SNK383_20_12810118_G_A_SNK383_1	20	12810118	G	A SNK383	145741	5200	1.6912E-05	4.04662E-05	1.6208E-09	4.02592E-05	0.000100009	0.0356797
185SNK384_20_12810118_G_A_SNK384_2	20	12810118	G	A SNK384	147601	5355	1.55507E-05	3.81814E-05	1.4429E-09	3.79861E-05	9.43628E-05	0.0362802
186SNK385_20_12810118_G_A_SNK385_3	20	12810118	G	A SNK385	144097	5336	2.44978E-05	8.29819E-05	6.815E-09	8.25531E-05	0.000205073	0.0370306
188SK215_5_73717969_G_A_SK215_1	5	73717969	G	A SK215	145975	16363	2.01463E-05	4.53848E-05	2.0383E-09	4.51478E-05	0.000112153	0.112095
205SNK312_5_173266954_G_A_SNK312_2	5	173266954	G	A SNK312	72517	1547	3.26982E-05	7.13994E-05	5.028E-09	7.09087E-05	0.000176147	0.0213329
17NA_5_174228431_G_C_S3PFC_1	5	174228431	G	C NA	49919	1534	2.48992E-05	4.02286E-05	1.6017E-09	4.2453E-05	0.000105459	0.015787505
18NA_5_174228431_G_C_S3PFC_3	5	174228431	G	C NA	34311	29	2.84456E-05	4.49986E-05	2.0029E-09	4.2453E-05	0.000105459	0.015787505
19NA_7_283913_T_A_S3PFC_2	7	283913	T	A NA	18171	0	5.29407E-05	0.000123851	1.5172E-08	0.000172125	0.000427581	0
24NA_9_136638046_C_T_S3PFC_1	9	136638046	C	T NA	47371	4	3.05441E-05	6.73903E-05	4.4756E-09	8.8965E-05	0.000221001	6.6676E-05
25NA_9_136638046_C_T_S3PFC_2	9	136638046	C	T NA	41069	2	3.5752E-05	6.04299E-05	3.5989E-09	8.8965E-05	0.000221001	6.6676E-05
26NA_9_136638046_C_T_S3PFC_3	9	136638046	C	T NA	29900	2	5.48824E-05	0.000126096	1.567E-08	8.8965E-05	0.000221001	6.6676E-05
74NA_2_17125698_C_T_S3PFC_1	2	17125698	C	T NA	28148	3	6.7018E-05	0.000117751	1.3709E-08	9.74024E-05	0.000241961	4.59349E-05
75NA_2_17125698_C_T_S3PFC_2	2	17125698	C	T NA	30646	0	4.8479E-05	9.22613E-05	8.4165E-09	9.74024E-05	0.000241961	4.59349E-05
76NA_2_17125698_C_T_S3PFC_3	2	17125698	C	T NA	32026	1	5.1095E-05	8.00484E-05	6.3357E-09	9.74024E-05	0.000241961	4.59349E-05
103NA_6_79286753_T_C_S3PFC_2	6	79286753	T	C NA	40164	0	4.95544E-05	9.39536E-05	8.727E-09	7.51117E-05	0.000186588	0
104NA_6_79286753_T_C_S3PFC_3	6	79286753	T	C NA	44679	0	3.4306E-05	5.0825E-05	2.5563E-09	7.51117E-05	0.000186588	0
111NA_8_40724674_G_A_S3PFC_1	8	40724674	G	A NA	24021	104	4.25707E-05	9.05005E-05	8.105E-09	0.000167721	0.000416643	0.004739067
112NA_8_40724674_G_A_S3PFC_2	8	40724674	G	A NA	21527	131	4.44042E-05	7.90849E-05	6.1899E-09	0.000167721	0.000416643	0.004739067
114NA_9_103459386_G_A_S3PFC_1	9	103459386	G	A NA	32109	83	5.84782E-05	7.6984E-05	5.8654E-09	6.59462E-05	0.00016382	0.00265107

TABLE 5-continued

Observed alternate allele fractions														
PrimerID	Chr	AlleleStart	Ref	Alt	Gene	IT Read Depth	Alt Allele Depth	Background AAF (within 50 nts)	Stdev AAF	Variance Background AAF	Stdev of Average Background		Confidence interval of Background	Average AAF
											Background	Background		
115NA_9_103459386_G_A_S3PFC_2	9	103459386	G	A	NA	30386	83	5.09081E-05	6.53611E-05	4.2285E-09	6.59462E-05	0.00016382	0.00265107	
116NA_9_103459386_G_A_S3PFC_3	9	103459386	G	A	NA	86091	227	6.13687E-05	5.46191E-05	2.9528E-09	6.59462E-05	0.00016382	0.00265107	
6_153444080_T_C_S3PFC_3	6	153444080	T	C	NA	58976	9758	0.00030193	8.73841E-05	7.5108E-09	7.83227E-05	0.000194564	0.44448	
5_174228431_G_C_S3PFC_1	5	174228431	G	C	NA	148622	149	1.24541E-05	2.27254E-05	5.1083E-10	2.26958E-05	5.63796E-05	0.000685806	
18NA_5_174228431_G_C_S3PFC_3	5	174228431	G	C	NA	75866	28	1.07793E-05	2.2912E-05	5.1937E-10	2.26958E-05	5.63796E-05	0.000685806	
7_283913_T_A_S3PFC_2	7	283913	T	A	NA	25754	0	1.86353E-05	4.34137E-05	1.864E-09	6.82206E-05	0.000169469	1.76498E-05	
7_283913_T_A_S3PFC_3	7	283913	T	A	NA	28329	1	3.50616E-05	8.67623E-05	7.4441E-09	6.82206E-05	0.000169469	1.76498E-05	
24NA_9_136638046_C_T_S3PFC_1	9	136638046	C	T	NA	86137	0	2.67017E-05	5.92953E-05	3.465E-09	7.96953E-05	0.000197974	4.74267E-05	
25NA_9_136638046_C_T_S3PFC_2	9	136638046	C	T	NA	135595	17	2.73487E-05	7.62347E-05	5.7275E-09	7.96953E-05	0.000197974	4.74267E-05	
26NA_9_136638046_C_T_S3PFC_3	9	136638046	C	T	NA	59147	1	3.56968E-05	0.000100033	9.8615E-09	7.96953E-05	0.000197974	4.74267E-05	
88NA_22_37475065_G_A_S3PFC_1	22	37475065	G	A	NA	24231	1	6.34762E-05	0.000128382	1.6171E-08	0.000121471	0.000301752	1.37565E-05	
89NA_22_37475065_G_A_S3PFC_2	22	37475065	G	A	NA	62505	0	2.41034E-05	5.70967E-05	3.225E-09	0.000121471	0.000301752	1.37565E-05	
90NA_22_37475065_G_A_S3PFC_3	22	37475065	G	A	NA	34638	0	7.69682E-05	0.000159271	2.487E-08	0.000121471	0.000301752	1.37565E-05	
103NA_6_79286753_T_C_S3PFC_2	6	79286753	T	C	NA	121355	0	4.01043E-05	7.43887E-05	5.4708E-09	5.96301E-05	0.000148129	0	
104NA_6_79286753_T_C_S3PFC_3	6	79286753	T	C	NA	86795	0	2.11582E-05	4.0716E-05	1.6407E-09	5.96301E-05	0.000148129	0	
111NA_8_40724674_G_A_S3PFC_1	8	40724674	G	A	NA	77763	153	2.25222E-05	5.78709E-05	3.3142E-09	6.513E-05	0.000161792	0.0019613	
112NA_8_40724674_G_A_S3PFC_2	8	40724674	G	A	NA	66376	108	2.61608E-05	7.10426E-05	4.995E-09	6.513E-05	0.000161792	0.0019613	
113NA_8_40724674_G_A_S3PFC_3	8	40724674	G	A	NA	111825	256	3.32318E-05	6.68025E-05	4.4166E-09	6.513E-05	0.000161792	0.0019613	
114NA_9_103459386_G_A_S3PFC_1	9	103459386	G	A	NA	43258	251	1.78852E-05	3.02772E-05	9.0726E-10	3.42549E-05	8.50939E-05	0.00484131	
115NA_9_103459386_G_A_S3PFC_2	9	103459386	G	A	NA	68024	264	1.84993E-05	3.55257E-05	1.2489E-09	3.42549E-05	8.50939E-05	0.00484131	
116NA_9_103459386_G_A_S3PFC_3	9	103459386	G	A	NA	76644	371	2.09892E-05	3.71224E-05	1.364E-09	3.42549E-05	8.50939E-05	0.00484131	
120NA_6_153444080_T_C_S3PFC_1	6	153444080	T	C	NA	38726	22027	1.20121E-05	2.65412E-05	6.9089E-10	7.80641E-05	0.000193922	0.30980405	
121NA_6_153444080_T_C_S3PFC_3	6	153444080	T	C	NA	105378	5355	2.47757E-05	0.000108161	1.1497E-08	7.80641E-05	0.000193922	0.30980405	

Example 2: Detecting Alleles with an Alternate Allele Fraction (AAF) at or Above 0.025%

[0143] To identify low frequency genetic variation in a target nucleic acid sequence with an alternate allele fraction of 0.025% or greater, three pairs of primers were designed to yield overlapping amplicons. Each pair of primers comprised a forward and a reverse primer, with each primer having a nucleotide sequence complementary to a portion of the target nucleic acid sequence. Each primer had an adapter at or near its 5' terminus and upstream from its complementary nucleic acid sequence. The adapter's nucleic acid sequence was complementary to a nucleic acid sequence used in an NGS platform, such as Ion Torrent or Illumina's MiSeq. Each individual reverse primer further comprised an index sequence upstream from the primer's complementary nucleic acid sequence. Additionally, each individual forward or reverse primer in each pair of primers further comprised a unique molecular identifier (UMI). No two primers had the same UMI.

[0144] Three distinct amplification reactions were prepared, each comprising one of the three pairs of primers. The reactions comprised 1.0 μ M primers, 1 $\times$ final concentration of 5 $\times$ Phusion High-fidelity Buffer (NEB), 200 μ M deoxynucleotide triphosphates (dNTPs), 0.1 μ l of 0.4 mM Biotin-14-dCTP, 1.0 units of Phusion High-fidelity Polymerase (NEB), and about 25 to 50 ng of template DNA. The reactions were subjected to an initial denaturation step of 30 seconds at 98° C. followed by 8 cycles of 98° C. (denaturing the template DNA) 10 seconds, 62° C. (annealing the primers to the template nucleic acid) for 20 seconds, and 72° C. (to extend the DNA product) for 30 seconds. After cycling, the reactions were subjected to an additional 10 minutes at 72° C. as a final extension step. The reaction products, or amplicons, were purified by washing 5 μ l of MyOne C1 streptavidin beads two times with 1 $\times$ Binding-Washing (B&W) buffer and then resuspending the beads in 25 μ l of 2 $\times$ B&W buffer. 25 μ l of the MyOne C1 streptavidin beads was then added to 25 μ l of the PCR amplicon and incubated at room temperature for 15 minutes with mixing. The mixture was exposed to a magnet, which isolates the beads with the amplicons bound thereto. The supernatant was removed, and 500 μ l 1 $\times$ B&W buffer was added to the beads, mixed, and exposed to the magnet. Again, the supernatant is removed, and the wash was repeated. The beads were finally resuspended in 28 μ l water. Some reaction products were purified using an exonuclease 1/shrimp alkaline phosphatase (ExoSap) enzymatic purification protocol, wherein 8 μ l of the commercially available ExoSap-It reagent (ThermoFisher) was added to the 20 μ l amplification reaction and incubated at 37° C. for 15 minutes followed by 80° C. for 15 minutes.

[0145] While the amplicons were attached to the streptavidin beads, an additional amplification was performed to enhance the copy number of the bound amplicons. Briefly, the additional amplification reactions comprised 1.0 μ M primers, 1 $\times$ final concentration of 5 $\times$ Phusion High-fidelity Buffer (NEB), 200 μ M deoxynucleotide triphosphates (dNTPs), 0.1 μ l of 0.4 mM Biotin-14-dCTP, 1.0 units of Phusion High-fidelity Polymerase (NEB), and about 25 to 50 ng of template DNA. The reactions were subjected to an initial denaturation step of 30 seconds at 98° C. followed by 20 cycles of 98° C. (denaturing the template DNA) 10 seconds, 62° C. (annealing the primers to the template nucleic acid) for 20 seconds, and 72° C. (to extend the DNA

product) for 30 seconds. After cycling, the reactions were subjected to an additional 10 minutes at 72° C. as a final extension step, and 5 μ l of the PCR reactions were pooled. A ThermoFisher MagJet purification kit that removes products <100 base pairs in length was used to purify the amplicons. Specifically, the amplicons in the pooled reactions were bound to streptavidin beads, and the supernatant was removed. The beads were then resuspended in 200 μ l of water, mixed, and incubated for two minutes. The mixture was then exposed to a magnet for two minutes, and the eluted DNA was captured.

[0146] Referring to FIGS. 11D to 11G, 1 μ l aliquots of eluted amplicons prepared using two rounds of amplification were run on a Bioanalyzer 2100 to confirm the quality of amplicons for use in downstream sequencing. FIG. 11D (first round=8 cycles; second round=20 cycles; biotin purification), FIG. 11E (first round=10 cycles; second round=20 cycles; biotin purification), FIG. 11F (first round=10 cycles; second round=20 cycles; no biotin purification), and FIG. 11G (first round=8 cycles; second round=25 cycles; ExoSAP purification) all show detectable amounts of the desired amplicons. For comparison purposes, data from an amplicon analyzed using TapeStation is shown in FIG. 12. Less sensitive than the Bioanalyzer 2100, the amplicons detected using the TapeStation are represented by much broader and rounded peaks compared to the Bioanalyzer 2100. However, this approach is still viable for the methods presented herein.

[0147] After determining the concentration of the eluted DNA, it was diluted to 100 pM, and the purified PCR reaction products were sequenced using the Ion Torrent system (ThermoFisher Scientific).

Example 3: Sensitivity and Reproducibility Assessment

[0148] The sensitivity and reproducibility of the methods described herein were assessed through serial dilutions of known germline mutations and known somatic mutations across a spectrum of alternative allele fractions. A comparison of alternative allele fractions with other known detections strategies including whole genome sequencing, whole exome sequencing, targeted sequencing, Sanger sequencing with Topo-cloning, and ddPCR was performed. First, triplicate primers (i.e., 3 unique pairs of primers) were designed as described in the methods for known germline mutations occurring in both the autosomal and X-chromosomal regions, including both heterozygous and hemizygous alleles. Twelve serial dilutions were sequenced on the Ion Torrent S5 with 400 base pair reads using six unique barcodes per primer. All reads were processed using custom analytical scripts (described in methods), allowing for the comparison of assessed and expected allelic fractions.

[0149] Referring to FIGS. 13 and 14, the methods described herein accurately measured alternative allele fractions as low as 0.025% and up to germline events when using a 50 ng of genomic DNA, although for significant detection above the amplicon-specific error rates, alternative allele fractions were typically required to be above 0.05%. The strong correlation between the expected and assessed alternative allele fractions ($R^2=0.9995$ and $R^2=0.9761$ for dilutions between 0-60% and for dilutions between 0-0.864%, respectively) across the assessed germline alleles, indicates that this method is extremely accurate for low-level alternative allele fractions.

[0150] Given that input DNA is often limited but is also known as an important factor for sensitivity for somatic alleles, decreased inputs of DNA were tested to determine if they could achieve a similar level of precision under the same dilution curve. Indeed, while decreased input DNA does impact the sensitivity, alternative allele fractions down to 0.05% remain detectable, though at a slightly elevated standard deviation among the triplicate primes for the lowest alternative allele fractions of 0.05%, indicating that when validating alleles below 0.1% alternative allele fractions, increased input DNA could improve precision. Furthermore, the impact of total sequencing depth on the accuracy was assessed to identify the minimum depth needed for accurate determination of alternative allele fractions. Using random sampling of the initial raw unmapped data, a strong correlation of read depths above threshold level can be made, and sequencing beyond this threshold will provide minimal benefits on the precision of the alternative allele fraction assessment.

Example 4: Somatic Mosaics in Human Brain Samples

[0151] Frozen postmortem human brain specimens from 61 autism spectrum disorder cases and 15 neurotypical controls were obtained for analysis. DNA was extracted from dorsolateral prefrontal cortex where available (or generic cortex in a minority of cases) using lysis buffer from the QIAamp DNA Mini kit (Qiagen) followed by phenol chloroform extraction and isopropanol cleanup. Samples UMB4334, UMB4899, UMB4999, UMB5027, UMB5115, UMB5176, UMB5297, UMB5302, UMB1638, UMB4671, and UMB797 were processed using TruSeq Nano DNA library preparation (Illumina) followed by Illumina HiSeq X Ten sequencing to a minimum 200× depth. All remaining samples were processed using TruSeq DNA PCR-Free library preparation (Illumina) followed by minimum 30× sequencing of seven separate libraries on the Illumina HiSeq X Ten, for a total minimum coverage of 210× per sample. An average of 251× depth was achieved across all samples, using 150 base pair paired-end reads. Two samples, UMB5771 and UMB5939, had parental saliva-derived DNA available, and DNA from both parents for these two cases was obtained and sequenced to about 50× depth. Parental DNA was not available for any other samples. Additionally, DNA was extracted from Brodmann Area 17 (occipital lobe) for cases UMB4638 and UMB4643 and sequenced at Macrogen to a minimum 210× depth following PCR-free library preparation. Bulk heart and liver sequencing data, as well as single-cell sequencing data from three individuals (UMB1465, UMB4643, and UMB4638) were used in this study.

Mutation Calling and Filtration

[0152] All paired-end FASTQ files were aligned using BWA-MEM version 0.7.8 to the GRCh37 human reference genome including the hs37d5 decoy sequence from the Broad Institute, following GATK best practices (software, broadinstitute.org/gatk/best-practices/). Mutect2-PoN was used to generate two pairs of panel-of-normals (PoN) by using 60 autism spectrum disorder samples or 15 control samples to remove sequencing artifacts and germline variants from the other group. Rare variants were further selected by filtering out any variant with a maximum popu-

lation minor allele frequency >0.001 in any of Kaviar, 1000 Genomes, EVS6500 (evs.gs.washington.edu/EVS/), ExACnonpsych, or gnomAD (gnomad.broadinstitute.org/). Repetitive region variants were removed using RepeatMasker (www.repeatmasker.org/), and variants within segmental duplication regions or shared between multiple individuals were also removed. Low-quality calls tagged “t_lod_fstar,” “str_contraction,” and “triallelic_site” were removed. For analysis of damaging heterozygous variants, variants were identified in the 78 risk genes previously used.

[0153] For somatic mutation detection, a minimum alternate (or variant) allele fraction (AAF or VAF) of 0.03 was required unless a variant was phasable by Mutect2, which allowed for rescue of variants down to an alternate allele fraction of 0.02. Low-quality calls tagged “triallelic_site” were removed. A minimum alternate read depth of four reads was required. Only private events among the population were analyzed. An upper alternate allele fraction threshold of 0.40 was set and heterozygous germline variants were removed. Variants within repetitive regions were also removed, leaving 14,984 candidate somatic mutations. MosaicForecast was then used to perform read-backed phasing and identify high-confidence mosaics from the candidate call set. Briefly, features likely to be correlated with mosaic detection specificity were selected: mapping quality, base quality, clustering of mutations, read depth, number of mismatches per read, read1/read2 bias, strand bias, base position, read position, trinucleotide context, sequencing cycle, library preparation method, and genotype likelihood. Based on these features a random forest model was trained using phased variants. Further training was conducted using parental whole genome sequencing data from two cases UMB5771 and UMB5939 as well as single cell whole genome sequencing data from three control brains, UMB1465, UMB4643, and UMB4638 for which inherited germline mutations or variants present in multiple single cells at a low alternate allele fraction (averaging alternate allele fraction <0.30, likely representing sequencing or alignment artifact), supplied a training set of false positives. Predicted mosaics were further filtered by removing genomic regions enriched for low-alternate allele fraction variants and by removing variants with unusually high sequencing depth that also occurred in regions marked as copy number variants (CNVs) by Meerkat. Following all training and filtration, 1143 putative mosaic variants were identified. One autism spectrum disorder sample, MSSM007, was eliminated from the study due to very high noise suggestive of contamination or sequencing artifact.

[0154] Pathogenicity prediction scores were calculated for functional mosaic and germline variants using SIFT, PolyPhen-2, MutationTaster, and CADD. To be considered damaging, a variant had to be predicted as damaging or probably damaging (or CADD phred score >20) by at least three out of four prediction tools. Mutations in genes were checked for overlap with the Simons Foundation Autism Research Initiative (SFARI) database of autism spectrum disorder—relevant genes (gene.sfari.org/), and with the Online Inheritance in Man (OMIM) database of genes with relevance to any human disease (www.omim.org/).

Triple Primer PCR Sequencing

[0155] Targeted validation was attempted on 243 of 1143 possible mosaic variants. PCR primers were designed for each variant and synthesized with Ion Torrent adapters P and

A, with barcodes added for unique identification. PCR amplification was performed using Phusion HotStart II DNA Polymerase (Thermo) as described by the manufacturer, with 20-25 cycles of amplification. Reactions were pooled and purified with AMPure XP technology (Agencourt), then sequenced on the Ion Torrent Personal Genome Machine using the Ion 530 chip with 400 base pair reads, reaching an average coverage of 118,000 reads per variant amongst reactions that yielded mappable reads. Following demultiplexing and trimming, reads were mapped using BWEMEM (a Burrows-Wheeler aligner algorithm) and locally realigned using GATK. BAM files were then imported into a CLC Genomics workbench (Qiagen) and mosaic variants were identified using the following filters: minimum frequency 0.05%, minimum depth 10,000× per reaction, minimum count 50, required significance 0.1%, central and neighborhood base quality of >15, and 3-nucleotide homopolymer filtration. Variants were then classified as validated true mosaics (198 variants), homozygous reference with variant not present (21 variants), germline heterozygous (1 variant), PCR reactions failed to amplify (19 variants), or undetermined (4 variants). The “undetermined” designation was used for variants for which the originally sequenced DNA was not available, so validation was conducted on a separate DNA extraction that could have slightly different clonal architecture. It was also used to classify two variants in which sequencing noise precluded validation interpretation. Validation success rates were calculated as the number of true mosaics divided by the sum of true mosaics, homozygous reference, and germline heterozygous. Weighted averaging across PCR and PCR-free variant validation was used to determine a comprehensive validation rate of 93%. Five variants from UMB5771 and UMB5939 were also re-sequenced in parent DNA, which confirmed a mosaic state in the offspring and homozygous reference in parents.

[0156] A deleterious missense C to A change in the autism spectrum disorder risk gene CACNA1A was called in 5.2% of sequencing reads in case UMB1174 (FIG. 15). Targeted validation of this region using the methods described herein generated 93,000 reads that confirmed an alternate allele fraction of 5.0%, meaning that this mutation is present in about 10% of cells.

[0157] Ion Torrent amplicon resequencing for 34 germline heterozygous mutations revealed that alternate allele frequencies were slightly over-dispersed compared to a binomial distribution (FIG. 16), likely due to noise induced by PCR amplification. The alternate allele frequency distribution was fit with a beta-binomial model to capture the over-dispersion ($\theta=452.44$, $p=1/(1+\theta)=0.0022$). 220 Ion Torrent-validated mosaics was used with a similar model to measure potential asymmetrical cell contributions to the brain during early embryonic development (FIG. 17A). Briefly, α_1 and $1-\alpha_1$ were defined as the fraction of brain cells deriving from each of the two cells created by the first division of the brain ancestor cell. A contribution parameter value of $\alpha=0.5$ meant that the first two cells contributed equally to the brain, while a non-0.5 value meant that the cell contribution was asymmetrical. Given a specific α_1 , it was possible to calculate the expected alternate allele frequency for mutations acquired at different branches of the early phylogeny (FIG. 17B). Assuming the mutation rate per cell generation was constant (i.e., the two cell divisions from the 2nd cell generation had the same mutation rate), the likeli-

hood of a mosaic arising on a specific branch was computed by multiplying the estimated sensitivity for detecting mosaics at the expected branch alternate allele frequency with the over-dispersion beta-binomial likelihood of the mosaic alternate allele fraction measured by the deep Ion Torrent sequencing. The log likelihoods for all sites were then summed over all branches to estimate the log likelihood of a specific al. al was fit by maximizing the log likelihood over $\alpha_1 \in [0.5, 1]$ using a grid search with step size=0.001. A likelihood ratio test was used to compare the asymmetrical model to the symmetrical model (i.e., $\alpha_1=0.5$), which clearly favored the model with unequal cell contribution during the 1st cell generation ($p<10^{-15}$). There is some evidence for asymmetrical contributions for later cell generations; however, since the asymmetric parameter α_1 estimated from the 2nd cell generation showed poor stability (FIG. 17C, $p=0.004$ compared to only one asymmetric cell division), asymmetric contribution was only assumed for the first cell generation. A 95% C.I. ([0.582, 0.607], FIG. 17D) was constructed using the likelihood ratio.

Example 5: Ultra-Sensitive Rapid Detection and Validation of Low-Frequency Somatic Mutations

[0158] The triple-primer PCR sequencing method substantially increases the throughput and sensitivity for the detection and validation of somatic mutations (FIGS. 4 and 5). This method utilizes multiple unique, carefully designed, custom primers targeting a region of interest in the genome to identify a novel mutation or assess the alternate allele fraction (AAF) of a known mutation in one or more samples. Unlike existing methods such as ddPCR, triple-primer PCR sequencing often requires little to no optimization after primer design and is less sensitive to DNA source, concentration, and nucleotide context. The robust sensitivity of the method detects and validates somatic and germline mutations using the Ion Torrent S5 platform and detects of novel alleles through modifications for Illumina sequencing.

Description of Triple Primer PCR Sequencing

[0159] While numerous studies have sought to define the error rates for the Ion Torrent platform due to the potential increased rate of insertion and deletion errors, particularly at homopolymers, the exact error rate appears to vary from sample to sample. Even more, while the rate of indel errors is likely elevated in the Ion Torrent platform over Illumina technology, the rates of SNV errors appear to be similar. It is likely that many estimates of errors are compounded by the combined effects of polymerase induced errors, mapping issues, and sequencing artifacts, all of which are known to reduce the sensitivity of detecting somatic mutations present in low fractions of a sample. Therefore, triple-primer PCR sequencing was developed to assess and partially mitigate these errors, while leveraging the rates to provide statistical confidence about a given mutation.

[0160] Prior studies have demonstrated the method of validating low AAF alleles using ultra-deep amplicon sequencing. However, technical issues including allelic dropout, artifacts (e.g., PCR- and sequencing platform-induced) and PCR duplicates can reduce the accuracy detected AAFs and possible result in both false negative calls as well as skewed AAFs. Triple-primer PCR sequencing overcomes these limitations through the use of multiple unique primers that are specifically designed to prevent

sharing binding sites while avoiding known mutations (i.e., individual specific and general population) but are within 250 nucleotides (nts) of the target mutation. Once designed, unique primer-specific barcodes are appended to the reverse primers, along with Ion Torrent adapters. Optionally, Illumina adapters and/or 10 nt molecular barcodes can be appended to the primers to improve sensitivity or usage on the Illumina platform. Customized primers amplify targets including the mutation or region of interest using reduced cycling and minimal amounts of DNA, and amplification products are sequenced on either the Ion Torrent S5 or Illumina MiSeq platform for ultra-deep coverage. This optimized process allows for independent analyses of each primer pair, determination of error rates bases on amplicon-specific error rates (i.e., level of PCR and sequencing induced artifacts across the amplicon), identification of allelic imbalances from additional mutations affecting primer binding or chromatin structure, and the assessment of the variation in AAF among primers. Together, these steps provide a robust and low-cost strategy for extremely precise estimation of AAFs which is broadly applicable to studies of somatic and germline mutations.

Accounting for Error Rates in Ion Torrent Data.

[0161] As the utility of the presently described invention relies on overcoming the previously described limitations of somatic mutation detection, triplicate unique primer sets were first designed around 5 known germline mutations (Tables 6A-6C) previously identified in bulk genomic DNA for testing the error rates of the method. The reduced PCR cycling conditions with a high-fidelity polymerase (4.4×10^{-7} ; Phusion HS, ThermoFisher) is estimated to result in an error rate of 8.8×10^{-6} at any given nucleotide position (ThermoFisher PCR Fidelity Calculator). Given that error

rates vary amongst amplicons due to the specific nucleotide content of each amplicon, an internal control was designed for assigning the significance of each identified mutation. Using these primers, background error rates from PCR and sequencing, the sensitivity to detect extremely low AAFs, accuracy of the ascertained AAF measurement, and required DNA input and sequencing depths were assessed.

[0162] First, reads and nucleotides were stringently filtered for nucleotide and mapping qualities ($q > 20$ and $Q > 20$), resulting in the removal of an average of 10% of bases at any given nucleotide position. Relaxing these parameters (e.g., $q10$, $Q10$) did not decrease the fraction of excluded sites or assessed AAF, supporting that most nucleotide positions are of high quality. Next, the rate of artifacts in the region of the amplicon surrounding the mutation of interest was assessed by the AAF of all alternate alleles at each position under the assumption that all non-reference high-quality alleles present at sites not known to have a mutation represent errors. Across all amplicons, a low average background mutation frequency ($0.018\% \text{ AAF} \pm 0.0067\%$) was found for nucleotides located in the flanking 50 nt on either site of a mutation. Consistent with prior studies, some amplicons exhibited positional variability in error rates due to mapping errors around indels, including artifacts arising during sequencing.

[0163] To further reduce the rate of indel-associated errors, a computational modeling approach that detects and corrects sequencing platform errors was incorporated. Specifically, Pollux, a recent error modeling algorithm that screens for and corrects an estimated $>95\%$ of all indel associated errors, was used. The correction of indel-associated errors resulted in nearly a 5-fold reduction in nucleotide error frequency ($0.0034\% \pm 0.0009\%$), allowing for mutations at extremely low AAFs to be distinguished from background sequencing and PCR-induced artifacts.

TABLE 6A

Chromosome	AlleleStart	AlleleEnd	Ref	Alt	Gene	Product Start	Product end	InsertStart	InsertEnd
X	153579431	153579431	T	C	FLNA	153579266	153579517	153579284	153579499
X	153579431	153579431	T	C	FLNA	153579289	153579555	153579311	153579536
X	153579431	153579431	T	C	FLNA	153579379	153579637	153579397	153579619
12	46321441	46321441	T	G	SCAF11	46321317	46321542	46321343	46321517
12	46321441	46321441	T	G	SCAF11	46321246	46321470	46321271	46321448
12	46321441	46321441	T	G	SCAF11	46321376	46321606	46321399	46321585
X	153594210	153594210	C	T	FLNA	153593965	153594295	153593983	153594277
X	153594210	153594210	C	T	FLNA	153594163	153594424	153594181	153594406
X	153594210	153594210	C	T	FLNA	153594114	153594378	153594132	153594360
16	3639306	3639306	G	A	SLX4	3639180	3639447	3639200	3639427
16	3639306	3639306	G	A	SLX4	3639109	3639337	3639129	3639319
16	3639306	3639306	G	A	SLX4	3639209	3639498	3639227	3639478
X	153599770	153599770	G	T	FLNA	153599611	153599868	153599629	153599850
X	153599770	153599770	G	T	FLNA	153599708	153599994	153599726	153599976
X	153599770	153599770	G	T	FLNA	153599747	153600008	153599766	153599989
18	21453038	21453038	C	T	LAMA3	21452938	21453163	21452959	21453143
18	21453038	21453038	C	T	LAMA3	21452848	21453097	21452867	21453076
18	21453038	21453038	C	T	LAMA3	21453007	21453231	21453025	21453208
X	153587777	153587777	G	C	FLNA	153587660	153587885	153587682	153587865
X	153587777	153587777	G	C	FLNA	153587508	153587801	153587528	153587781
X	153587777	153587777	G	C	FLNA	153587606	153587897	153587626	153587878
19	44571260	44571260	C	A	ZNF223	44571155	44571379	44571175	44571359
19	44571260	44571260	C	A	ZNF223	44571066	44571291	44571085	44571270
19	44571260	44571260	C	A	ZNF223	44571227	44571456	44571251	44571429
X	153579431	153579431	T	C	FLNA	153579266	153579517	153579284	153579499
X	153579431	153579431	T	C	FLNA	153579289	153579555	153579311	153579536
X	153579431	153579431	T	C	FLNA	153579379	153579637	153579397	153579619
12	46321441	46321441	T	G	SCAF11	46321317	46321542	46321343	46321517
12	46321441	46321441	T	G	SCAF11	46321246	46321470	46321271	46321448
12	46321441	46321441	T	G	SCAF11	46321376	46321606	46321399	46321585

TABLE 6A-continued

Chromosome	AlleleStart	AlleleEnd	Ref	Alt	Gene	Product Start	Product end	InsertStart	InsertEnd
X	153594210	153594210	C	T	FLNA	153593965	153594295	153593983	153594277
X	153594210	153594210	C	T	FLNA	153594163	153594424	153594181	153594406
X	153594210	153594210	C	T	FLNA	153594114	153594378	153594132	153594360
16	3639306	3639306	G	A	SLX4	3639180	3639447	3639200	3639427
16	3639306	3639306	G	A	SLX4	3639109	3639337	3639129	3639319
16	3639306	3639306	G	A	SLX4	3639209	3639498	3639227	3639478
X	153599770	153599770	G	T	FLNA	153599611	153599868	153599629	153599850
X	153599770	153599770	G	T	FLNA	153599708	153599994	153599726	153599976
X	153599770	153599770	G	T	FLNA	153599747	153600008	153599766	153599989
18	21453038	21453038	C	T	LAMA3	21452938	21453163	21452959	21453143
18	21453038	21453038	C	T	LAMA3	21452848	21453097	21452867	21453076
18	21453038	21453038	C	T	LAMA3	21453007	21453231	21453025	21453208
X	153587777	153587777	G	C	FLNA	153587660	153587885	153587682	153587865
X	153587777	153587777	G	C	FLNA	153587508	153587801	153587528	153587781
X	153587777	153587777	G	C	FLNA	153587606	153587897	153587626	153587878
19	44571260	44571260	C	A	ZNF223	44571155	44571379	44571175	44571359
19	44571260	44571260	C	A	ZNF223	44571066	44571291	44571085	44571270
19	44571260	44571260	C	A	ZNF223	44571227	44571456	44571251	44571429
X	153579431	153579431	T	C	FLNA	153579266	153579517	153579284	153579499
X	153579431	153579431	T	C	FLNA	153579289	153579555	153579311	153579536
X	153579431	153579431	T	C	FLNA	153579379	153579637	153579397	153579619
12	46321441	46321441	T	G	SCAF11	46321317	46321542	46321343	46321517
12	46321441	46321441	T	G	SCAF11	46321246	46321470	46321271	46321448
12	46321441	46321441	T	G	SCAF11	46321376	46321606	46321399	46321585
X	153594210	153594210	C	T	FLNA	153593965	153594295	153593983	153594277
X	153594210	153594210	C	T	FLNA	153594163	153594424	153594181	153594406
X	153594210	153594210	C	T	FLNA	153594114	153594378	153594132	153594360
16	3639306	3639306	G	A	SLX4	3639180	3639447	3639200	3639427
16	3639306	3639306	G	A	SLX4	3639109	3639337	3639129	3639319
16	3639306	3639306	G	A	SLX4	3639209	3639498	3639227	3639478
X	153599770	153599770	G	T	FLNA	153599611	153599868	153599629	153599850
X	153599770	153599770	G	T	FLNA	153599708	153599994	153599726	153599976
X	153599770	153599770	G	T	FLNA	153599747	153600008	153599766	153599989
18	21453038	21453038	C	T	LAMA3	21452938	21453163	21452959	21453143
18	21453038	21453038	C	T	LAMA3	21452848	21453097	21452867	21453076
18	21453038	21453038	C	T	LAMA3	21453007	21453231	21453025	21453208
X	153587777	153587777	G	C	FLNA	153587660	153587885	153587682	153587865
X	153587777	153587777	G	C	FLNA	153587508	153587801	153587528	153587781
X	153587777	153587777	G	C	FLNA	153587606	153587897	153587626	153587878
19	44571260	44571260	C	A	ZNF223	44571155	44571379	44571175	44571359
19	44571260	44571260	C	A	ZNF223	44571066	44571291	44571085	44571270
19	44571260	44571260	C	A	ZNF223	44571227	44571456	44571251	44571429
X	153579431	153579431	T	C	FLNA	153579266	153579517	153579284	153579499
X	153579431	153579431	T	C	FLNA	153579289	153579555	153579311	153579536
X	153579431	153579431	T	C	FLNA	153579379	153579637	153579397	153579619
12	46321441	46321441	T	G	SCAF11	46321317	46321542	46321343	46321517
12	46321441	46321441	T	G	SCAF11	46321246	46321470	46321271	46321448
12	46321441	46321441	T	G	SCAF11	46321376	46321606	46321399	46321585
X	153594210	153594210	C	T	FLNA	153593965	153594295	153593983	153594277
X	153594210	153594210	C	T	FLNA	153594163	153594424	153594181	153594406
X	153594210	153594210	C	T	FLNA	153594114	153594378	153594132	153594360
16	3639306	3639306	G	A	SLX4	3639180	3639447	3639200	3639427
16	3639306	3639306	G	A	SLX4	3639109	3639337	3639129	3639319
16	3639306	3639306	G	A	SLX4	3639209	3639498	3639227	3639478
X	153599770	153599770	G	T	FLNA	153599611	153599868	153599629	153599850
X	153599770	153599770	G	T	FLNA	153599708	153599994	153599726	153599976
X	153599770	153599770	G	T	FLNA	153599747	153600008	153599766	153599989
18	21453038	21453038	C	T	LAMA3	21452938	21453163	21452959	21453143
18	21453038	21453038	C	T	LAMA3	21452848	21453097	21452867	21453076
18	21453038	21453038	C	T	LAMA3	21453007	21453231	21453025	21453208
X	153587777	153587777	G	C	FLNA	153587660	153587885	153587682	153587865
X	153587777	153587777	G	C	FLNA	153587508	153587801	153587528	153587781
X	153587777	153587777	G	C	FLNA	153587606	153587897	153587626	153587878
19	44571260	44571260	C	A	ZNF223	44571155	44571379	44571175	44571359
19	44571260	44571260	C	A	ZNF223	44571066	44571291	44571085	44571270
19	44571260	44571260	C	A	ZNF223	44571227	44571456	44571251	44571429
X	153579431	153579431	T	C	FLNA	153579266	153579517	153579284	153579499
X	153579431	153579431	T	C	FLNA	153579289	153579555	153579311	153579536
X	153579431	153579431	T	C	FLNA	153579379	153579637	153579397	153579619
12	46321441	46321441	T	G	SCAF11	46321317	46321542	46321343	46321517
12	46321441	46321441	T	G	SCAF11	46321246	46321470	46321271	46321448
12	46321441	46321441	T	G	SCAF11	46321376	46321606	46321399	46321585
X	153594210	153594210	C	T	FLNA	153593965	153594295	153593983	153594277
X	153594210	153594210	C	T	FLNA	153594163	153594424	153594181	153594406

TABLE 6A-continued

Chromosome	AlleleStart	AlleleEnd	Ref	Alt	Gene	Product Start	Product end	InsertStart	InsertEnd
X	153594210	153594210	C	T	FLNA	153594114	153594378	153594132	153594360
16	3639306	3639306	G	A	SLX4	3639180	3639447	3639200	3639427
16	3639306	3639306	G	A	SLX4	3639109	3639337	3639129	3639319
16	3639306	3639306	G	A	SLX4	3639209	3639498	3639227	3639478
X	153599770	153599770	G	T	FLNA	153599611	153599868	153599629	153599850
X	153599770	153599770	G	T	FLNA	153599708	153599994	153599726	153599976
X	153599770	153599770	G	T	FLNA	153599747	153600008	153599766	153599989
18	21453038	21453038	C	T	LAMA3	21452938	21453163	21452959	21453143
18	21453038	21453038	C	T	LAMA3	21452848	21453097	21452867	21453076
18	21453038	21453038	C	T	LAMA3	21453007	21453231	21453025	21453208
X	153587777	153587777	G	C	FLNA	153587660	153587885	153587682	153587865
X	153587777	153587777	G	C	FLNA	153587508	153587801	153587528	153587781
X	153587777	153587777	G	C	FLNA	153587606	153587897	153587626	153587878
19	44571260	44571260	C	A	ZNF223	44571155	44571379	44571175	44571359
19	44571260	44571260	C	A	ZNF223	44571066	44571291	44571085	44571270
19	44571260	44571260	C	A	ZNF223	44571227	44571456	44571251	44571429
X	153579431	153579431	T	C	FLNA	153579266	153579517	153579284	153579499
X	153579431	153579431	T	C	FLNA	153579289	153579555	153579311	153579536
X	153579431	153579431	T	C	FLNA	153579379	153579637	153579397	153579619
12	46321441	46321441	T	G	SCAF11	46321317	46321542	46321343	46321517
12	46321441	46321441	T	G	SCAF11	46321246	46321470	46321271	46321448
12	46321441	46321441	T	G	SCAF11	46321376	46321606	46321399	46321585
X	153594210	153594210	C	T	FLNA	153593965	153594295	153593983	153594277
X	153594210	153594210	C	T	FLNA	153594163	153594424	153594181	153594406
X	153594210	153594210	C	T	FLNA	153594114	153594378	153594132	153594360
16	3639306	3639306	G	A	SLX4	3639180	3639447	3639200	3639427
16	3639306	3639306	G	A	SLX4	3639109	3639337	3639129	3639319
16	3639306	3639306	G	A	SLX4	3639209	3639498	3639227	3639478
X	153599770	153599770	G	T	FLNA	153599611	153599868	153599629	153599850
X	153599770	153599770	G	T	FLNA	153599708	153599994	153599726	153599976
X	153599770	153599770	G	T	FLNA	153599747	153600008	153599766	153599989
18	21453038	21453038	C	T	LAMA3	21452938	21453163	21452959	21453143
18	21453038	21453038	C	T	LAMA3	21452848	21453097	21452867	21453076
18	21453038	21453038	C	T	LAMA3	21453007	21453231	21453025	21453208
X	153587777	153587777	G	C	FLNA	153587660	153587885	153587682	153587865
X	153587777	153587777	G	C	FLNA	153587508	153587801	153587528	153587781
X	153587777	153587777	G	C	FLNA	153587606	153587897	153587626	153587878
19	44571260	44571260	C	A	ZNF223	44571155	44571379	44571175	44571359
19	44571260	44571260	C	A	ZNF223	44571066	44571291	44571085	44571270
19	44571260	44571260	C	A	ZNF223	44571227	44571456	44571251	44571429

TABLE 6B

Chromosome	AlleleStart	AlleleEnd	Forward
X	153579431	153579431	CCTCTCTATGGGCAGTCGGTGATCAGGGCCTCACCTTGCTC
X	153579431	153579431	CCTCTCTATGGGCAGTCGGTGATCTGTGACATAGCACTCCTCCAG
X	153579431	153579431	CCTCTCTATGGGCAGTCGGTGATAGGCTGGCTGGTTGACCT
12	46321441	46321441	CCTCTCTATGGGCAGTCGGTGATAATCACACTCCATAGGTATCATTTC
12	46321441	46321441	CCTCTCTATGGGCAGTCGGTGATTTCATTTCATTGTTTAAAGATCAGCA
12	46321441	46321441	CCTCTCTATGGGCAGTCGGTGATTCAATGTGTGTTTTAGGCAACTC
X	153594210	153594210	CCTCTCTATGGGCAGTCGGTGATAGGGGACATGCAAGACA
X	153594210	153594210	CCTCTCTATGGGCAGTCGGTGATGCGAGCTCTCCGAAGGT
X	153594210	153594210	CCTCTCTATGGGCAGTCGGTGATGTTGACCCTGTGGGCAGA
16	3639306	3639306	CCTCTCTATGGGCAGTCGGTGATTCCCTTGCGTAGTGCAGCTT
16	3639306	3639306	CCTCTCTATGGGCAGTCGGTGATCAGAGCCGAATCAGAAAGC
16	3639306	3639306	CCTCTCTATGGGCAGTCGGTGATGGGGTGGTGTCAGAGAT
X	153599770	153599770	CCTCTCTATGGGCAGTCGGTGATCATTTTGAGGCGCGAGAA

TABLE 6B-continued

Chromosome	AlleleStart	AlleleEnd	Forward
X	153599770	153599770	CCTCTCTATGGGCAGTCGGTGATGAGGCAGGGAGCAGAGGT
X	153599770	153599770	CCTCTCTATGGGCAGTCGGTGATCCTTTAAATGCGGGAGGAG
18	21453038	21453038	CCTCTCTATGGGCAGTCGGTGATGAGCAGGAAGGGCAGGTATAA
18	21453038	21453038	CCTCTCTATGGGCAGTCGGTGATCTGGCACAGGCTGACTCAT
18	21453038	21453038	CCTCTCTATGGGCAGTCGGTGATGGATGCCCTCCAGCAGTGA
X	153587777	153587777	CCTCTCTATGGGCAGTCGGTGATAGCCTCATAAGGGATGTACTCG
X	153587777	153587777	CCTCTCTATGGGCAGTCGGTGATCCTTGAAAGGACTGCCTGAG
X	153587777	153587777	CCTCTCTATGGGCAGTCGGTGATCTCCTCACCTGGCACTTGAT
19	44571260	44571260	CCTCTCTATGGGCAGTCGGTGATAGAGCCACACAGGAGAGAG
19	44571260	44571260	CCTCTCTATGGGCAGTCGGTGATATCAGCGAGTCCACACTGG
19	44571260	44571260	CCTCTCTATGGGCAGTCGGTGATTGTAATCATAAGAGACTCCATTGC
X	153579431	153579431	CCTCTCTATGGGCAGTCGGTGATCAGGGCCTCACCTTGGTC
X	153579431	153579431	CCTCTCTATGGGCAGTCGGTGATCTGTGACATAGCACTCCTCCAG
X	153579431	153579431	CCTCTCTATGGGCAGTCGGTGATAGGCTGGCTGGTTGACCT
12	46321441	46321441	CCTCTCTATGGGCAGTCGGTGATAATCACACTCCATAGGTATCATTTCA
12	46321441	46321441	CCTCTCTATGGGCAGTCGGTGATTTCATTCATTTGTTTAAGATCAGCA
12	46321441	46321441	CCTCTCTATGGGCAGTCGGTGATTCAATGTGTGTTTATAGCAACTC
X	153594210	153594210	CCTCTCTATGGGCAGTCGGTGATAGGGGGACATGCAAGACA
X	153594210	153594210	CCTCTCTATGGGCAGTCGGTGATGCGAGCTCTTCCGAAGGT
X	153594210	153594210	CCTCTCTATGGGCAGTCGGTGATGTTGACCTGTGGGCAGA
16	3639306	3639306	CCTCTCTATGGGCAGTCGGTGATTCTCTGGGTAGTGACGCTT
16	3639306	3639306	CCTCTCTATGGGCAGTCGGTGATCAGAGCCGAATTGAGAAAGC
16	3639306	3639306	CCTCTCTATGGGCAGTCGGTGATGGGGTGGTGTCCAGGAGT
X	153599770	153599770	CCTCTCTATGGGCAGTCGGTGATCATTTTGGAGCGCGAGAA
X	153599770	153599770	CCTCTCTATGGGCAGTCGGTGATGAGGCAGGGAGCAGAGGT
X	153599770	153599770	CCTCTCTATGGGCAGTCGGTGATCCTTTAAATGCGGGAGGAG
18	21453038	21453038	CCTCTCTATGGGCAGTCGGTGATGAGCAGGAAGGGCAGGTATAA
18	21453038	21453038	CCTCTCTATGGGCAGTCGGTGATCTGGCACAGGCTGACTCAT
18	21453038	21453038	CCTCTCTATGGGCAGTCGGTGATGGATGCCCTCCAGCAGTGA
X	153587777	153587777	CCTCTCTATGGGCAGTCGGTGATAGCCTCATAAGGGATGTACTCG
X	153587777	153587777	CCTCTCTATGGGCAGTCGGTGATCCTTGAAAGGACTGCCTGAG
X	153587777	153587777	CCTCTCTATGGGCAGTCGGTGATCTCCTCACCTGGCACTTGAT
19	44571260	44571260	CCTCTCTATGGGCAGTCGGTGATAGAGCCACACAGGAGAGAG
19	44571260	44571260	CCTCTCTATGGGCAGTCGGTGATATCAGCGAGTCCACACTGG
19	44571260	44571260	CCTCTCTATGGGCAGTCGGTGATTGTAATCATAAGAGACTCCATTGC
X	153579431	153579431	CCTCTCTATGGGCAGTCGGTGATCAGGGCCTCACCTTGGTC
X	153579431	153579431	CCTCTCTATGGGCAGTCGGTGATCTGTGACATAGCACTCCTCCAG
X	153579431	153579431	CCTCTCTATGGGCAGTCGGTGATAGGCTGGCTGGTTGACCT

TABLE 6B-continued

Chromosome	AlleleStart	AlleleEnd	Forward
12	46321441	46321441	CCTCTCTATGGGCAGTCGGTGATAATCACACTCCATAGGTATCATTCA
12	46321441	46321441	CCTCTCTATGGGCAGTCGGTGATTTTCATTTCATTTGTTTAAGATCAGCA
12	46321441	46321441	CCTCTCTATGGGCAGTCGGTGATTCAATGTGTGTTTATAGGCAACTC
X	153594210	153594210	CCTCTCTATGGGCAGTCGGTGATAGGGGGACATGCAAGACA
X	153594210	153594210	CCTCTCTATGGGCAGTCGGTGATGCGAGCTCTTCCGAAGGT
X	153594210	153594210	CCTCTCTATGGGCAGTCGGTGATGTTGACCTGTGGGCAGA
16	3639306	3639306	CCTCTCTATGGGCAGTCGGTGATTCCTCTGGGTAGTGACGCTT
16	3639306	3639306	CCTCTCTATGGGCAGTCGGTGATCAGAGCCGAATTCAGAAAGC
16	3639306	3639306	CCTCTCTATGGGCAGTCGGTGATGGGGTGGTGTCCAGGAGT
X	153599770	153599770	CCTCTCTATGGGCAGTCGGTGATCATTTTGAGGCGCAGAGAA
X	153599770	153599770	CCTCTCTATGGGCAGTCGGTGATGAGGCAGGGAGCAGAGGT
X	153599770	153599770	CCTCTCTATGGGCAGTCGGTGATCCTTTAAATGCGGGAGGAG
18	21453038	21453038	CCTCTCTATGGGCAGTCGGTGATGAGCAGGAAGGGCAGGTATAA
18	21453038	21453038	CCTCTCTATGGGCAGTCGGTGATCTGGCACAGGCTGACTCAT
18	21453038	21453038	CCTCTCTATGGGCAGTCGGTGATGGATGCCTCCAGCAGTGA
X	153587777	153587777	CCTCTCTATGGGCAGTCGGTGATAGCCTCATAAGGGATGTACTCG
X	153587777	153587777	CCTCTCTATGGGCAGTCGGTGATCCTTGAAAGGACTGCCTGAG
X	153587777	153587777	CCTCTCTATGGGCAGTCGGTGATCTCCTCACCTGGCACTTGAT
19	44571260	44571260	CCTCTCTATGGGCAGTCGGTGATAGAGCCACACAGGAGAGAG
19	44571260	44571260	CCTCTCTATGGGCAGTCGGTGATATCAGCGAGTCCACACTGG
19	44571260	44571260	CCTCTCTATGGGCAGTCGGTGATTTGAATCATAAGAGACTCCATTGC
X	153579431	153579431	CCTCTCTATGGGCAGTCGGTGATCAGGGCCTCACCTTGGTGTC
X	153579431	153579431	CCTCTCTATGGGCAGTCGGTGATCTGTGACATAGCACTCCTCCAG
X	153579431	153579431	CCTCTCTATGGGCAGTCGGTGATAGGCTGGCTGGTTGACCT
12	46321441	46321441	CCTCTCTATGGGCAGTCGGTGATAATCACACTCCATAGGTATCATTCA
12	46321441	46321441	CCTCTCTATGGGCAGTCGGTGATTTTCATTTCATTTGTTTAAGATCAGCA
12	46321441	46321441	CCTCTCTATGGGCAGTCGGTGATTCAATGTGTGTTTATAGGCAACTC
X	153594210	153594210	CCTCTCTATGGGCAGTCGGTGATAGGGGGACATGCAAGACA
X	153594210	153594210	CCTCTCTATGGGCAGTCGGTGATGCGAGCTCTTCCGAAGGT
X	153594210	153594210	CCTCTCTATGGGCAGTCGGTGATGTTGACCTGTGGGCAGA
16	3639306	3639306	CCTCTCTATGGGCAGTCGGTGATTCCTCTGGGTAGTGACGCTT
16	3639306	3639306	CCTCTCTATGGGCAGTCGGTGATCAGAGCCGAATTCAGAAAGC
16	3639306	3639306	CCTCTCTATGGGCAGTCGGTGATGGGGTGGTGTCCAGGAGT
X	153599770	153599770	CCTCTCTATGGGCAGTCGGTGATCATTTTGAGGCGCAGAGAA
X	153599770	153599770	CCTCTCTATGGGCAGTCGGTGATGAGGCAGGGAGCAGAGGT
X	153599770	153599770	CCTCTCTATGGGCAGTCGGTGATCCTTTAAATGCGGGAGGAG
18	21453038	21453038	CCTCTCTATGGGCAGTCGGTGATGAGCAGGAAGGGCAGGTATAA

TABLE 6B-continued

Chromosome	AlleleStart	AlleleEnd	Forward
18	21453038	21453038	CCTCTCTATGGGCAGTCGGTGATCTGGCACAGGCTGACTCAT
18	21453038	21453038	CCTCTCTATGGGCAGTCGGTGATGGATGCCCTCCAGCAGTGA
X	153587777	153587777	CCTCTCTATGGGCAGTCGGTGATAGCCTCATAAGGGATGTACTCG
X	153587777	153587777	CCTCTCTATGGGCAGTCGGTGATCCTTGAAAGGACTGCCTGAG
X	153587777	153587777	CCTCTCTATGGGCAGTCGGTGATCTCCTCACCTGGCACTTGAT
19	44571260	44571260	CCTCTCTATGGGCAGTCGGTGATAGAGCCACACAGGAGAGAG
19	44571260	44571260	CCTCTCTATGGGCAGTCGGTGATATCAGCGAGTCCACACTGG
19	44571260	44571260	CCTCTCTATGGGCAGTCGGTGATTGAATCATAAGAGACTCCATTGC
X	153579431	153579431	CCTCTCTATGGGCAGTCGGTGATCAGGGCCTCACCTTGGTC
X	153579431	153579431	CCTCTCTATGGGCAGTCGGTGATCTGTGACATAGCACTCCTCCAG
X	153579431	153579431	CCTCTCTATGGGCAGTCGGTGATAGGCTGGCTGGTTGACCT
12	46321441	46321441	CCTCTCTATGGGCAGTCGGTGATAATCACACTCCATAGGTATCATTCA
12	46321441	46321441	CCTCTCTATGGGCAGTCGGTGATTTCATTTCATTGTTTAAGATCAGCA
12	46321441	46321441	CCTCTCTATGGGCAGTCGGTGATTCAATGTGTGTTTTAGGCAACTC
X	153594210	153594210	CCTCTCTATGGGCAGTCGGTGATAGGGGGACATGCAAGACA
X	153594210	153594210	CCTCTCTATGGGCAGTCGGTGATGCGAGCTCTCCGAAGGT
X	153594210	153594210	CCTCTCTATGGGCAGTCGGTGATGTTGACCCTGTGGGCAGA
16	3639306	3639306	CCTCTCTATGGGCAGTCGGTGATCCTCTGGGTAGTGCAGCTT
16	3639306	3639306	CCTCTCTATGGGCAGTCGGTGATCAGAGCCGAATTAGAAAGC
16	3639306	3639306	CCTCTCTATGGGCAGTCGGTGATGGGGTGGTGTCAGGAGT
X	153599770	153599770	CCTCTCTATGGGCAGTCGGTGATCATTTTGAGGCGCGAGAA
X	153599770	153599770	CCTCTCTATGGGCAGTCGGTGATGAGGCAGGGAGCAGAGGT
X	153599770	153599770	CCTCTCTATGGGCAGTCGGTGATCCTTTAAATGCGGGAGGAG
18	21453038	21453038	CCTCTCTATGGGCAGTCGGTGATGAGCAGGAAGGGCAGGTATAA
18	21453038	21453038	CCTCTCTATGGGCAGTCGGTGATCTGGCACAGGCTGACTCAT
18	21453038	21453038	CCTCTCTATGGGCAGTCGGTGATGGATGCCCTCCAGCAGTGA
X	153587777	153587777	CCTCTCTATGGGCAGTCGGTGATAGCCTCATAAGGGATGTACTCG
X	153587777	153587777	CCTCTCTATGGGCAGTCGGTGATCCTTGAAAGGACTGCCTGAG
X	153587777	153587777	CCTCTCTATGGGCAGTCGGTGATCTCCTCACCTGGCACTTGAT
19	44571260	44571260	CCTCTCTATGGGCAGTCGGTGATAGAGCCACACAGGAGAGAG
19	44571260	44571260	CCTCTCTATGGGCAGTCGGTGATATCAGCGAGTCCACACTGG
19	44571260	44571260	CCTCTCTATGGGCAGTCGGTGATTGAATCATAAGAGACTCCATTGC
X	153579431	153579431	CCTCTCTATGGGCAGTCGGTGATCAGGGCCTCACCTTGGTC
X	153579431	153579431	CCTCTCTATGGGCAGTCGGTGATCTGTGACATAGCACTCCTCCAG
X	153579431	153579431	CCTCTCTATGGGCAGTCGGTGATAGGCTGGCTGGTTGACCT
12	46321441	46321441	CCTCTCTATGGGCAGTCGGTGATAATCACACTCCATAGGTATCATTCA
12	46321441	46321441	CCTCTCTATGGGCAGTCGGTGATTTCATTTCATTGTTTAAGATCAGCA
12	46321441	46321441	CCTCTCTATGGGCAGTCGGTGATTCAATGTGTGTTTTAGGCAACTC

TABLE 6B-continued

Chromosome	AlleleStart	AlleleEnd	Forward
X	153594210	153594210	CCTCTCTATGGGCAGTCGGTGATAGGGGGACATGCAAGACA
X	153594210	153594210	CCTCTCTATGGGCAGTCGGTGATGCGAGCTCTTCCGAAGGT
X	153594210	153594210	CCTCTCTATGGGCAGTCGGTGATGTTGACCTGTGGGCAGA
16	3639306	3639306	CCTCTCTATGGGCAGTCGGTGATTCTCTGGGTAGTGCAGCTT
16	3639306	3639306	CCTCTCTATGGGCAGTCGGTGATCAGAGCCGAATTCAGAAAGC
16	3639306	3639306	CCTCTCTATGGGCAGTCGGTGATGGGGTGGTGTCCAGGAGT
X	153599770	153599770	CCTCTCTATGGGCAGTCGGTGATCATTTTGAGGCGCGAGAA
X	153599770	153599770	CCTCTCTATGGGCAGTCGGTGATGAGGCAGGGAGCAGAGGT
X	153599770	153599770	CCTCTCTATGGGCAGTCGGTGATCCTTTAAATGCGGGAGGAG
18	21453038	21453038	CCTCTCTATGGGCAGTCGGTGATGAGCAGGAAGGGCAGGTATAA
18	21453038	21453038	CCTCTCTATGGGCAGTCGGTGATCTGGCACAGGCTGACTCAT
18	21453038	21453038	CCTCTCTATGGGCAGTCGGTGATGGATGCCCTCCAGCAGTGA
X	153587777	153587777	CCTCTCTATGGGCAGTCGGTGATAGCCTCATAAGGGATGTACTCG
X	153587777	153587777	CCTCTCTATGGGCAGTCGGTGATCCTTGAAAGGACTGCCCTGAG
X	153587777	153587777	CCTCTCTATGGGCAGTCGGTGATCTCCTCACCTGGCACTTGAT
19	44571260	44571260	CCTCTCTATGGGCAGTCGGTGATAGAGCCACACAGGAGAGAG
19	44571260	44571260	CCTCTCTATGGGCAGTCGGTGATATCAGCGAGTCCACACTGG
19	44571260	44571260	CCTCTCTATGGGCAGTCGGTGATTGAATCATAAGAGACTCCATTGC

TABLE 6C

Chromosome	AlleleStart	AlleleEnd	Reverse	Barcode
X	153579431	153579431	CCATCTCATCCCTCGGTGTCCTCGACTCAGTtaacggacgCGCAGATGGTAAGTGC	ttaacggacg
X	153579431	153579431	CCATCTCATCCCTCGGTGTCCTCGACTCAGTcaggcttaeTGCAATCAGTGGCTCTCC	tccggcttac
X	153579431	153579431	CCATCTCATCCCTCGGTGTCCTCGACTCAGTcgtctcattcagTCCCTTCCTGCGCACCTG	tctcattcag
12	46321441	46321441	CCATCTCATCCCTCGGTGTCCTCGACTCAGGcggtcatacACATGTGATACTTTTGGGAATGAA G	gcggtcatac
12	46321441	46321441	CCATCTCATCCCTCGGTGTCCTCGACTCAGTaggacgttcCTTCTGAACACCAAAATTGGAAA	taggacgttc
12	46321441	46321441	CCATCTCATCCCTCGGTGTCCTCGACTCAGacgacgcaacTGTTAAGAGCCCGAGAGTTCA	acgacgcaac
X	153594210	153594210	CCATCTCATCCCTCGGTGTCCTCGACTCAGcttctcggacGGGGCCCTACTCTTTGA	cttctcggac
X	153594210	153594210	CCATCTCATCCCTCGGTGTCCTCGACTCAGctatcgcogttCTCGACGCCCTACACTG	cattgcggtta
X	153594210	153594210	CCATCTCATCCCTCGGTGTCCTCGACTCAGcagacagaaTGAAGTGCCTCTCTGCTG	tgaggacggc
16	3639306	3639306	CCATCTCATCCCTCGGTGTCCTCGACTCAGtgaggacggcAGTGACGATGAGCAGAGGT	gcctgcgcag
16	3639306	3639306	CCATCTCATCCCTCGGTGTCCTCGACTCAGgctgcgcagGCCAATTCOCATTGACCA	gttgacgtct
16	3639306	3639306	CCATCTCATCCCTCGGTGTCCTCGACTCAGgttgaegtctCCAAGTTTCGTGAACGAGAC	gagatcgatt
X	153599770	153599770	CCATCTCATCCCTCGGTGTCCTCGACTCAGctcagctcCTAGTGGGGGCATTCCAA	agttcgagcc
X	153599770	153599770	CCATCTCATCCCTCGGTGTCCTCGACTCAGagttcgagccCTTAGGGCGCGTTTCCT	ctcaggctca
X	153599770	153599770	CCATCTCATCCCTCGGTGTCCTCGACTCAGggcaataataTCCACATACCTCGCTTGCA	ggcaataata
18	21453038	21453038	CCATCTCATCCCTCGGTGTCCTCGACTCAGggtactcctgCAACTGTAGCCCGACACTGC	ggtagctcag
18	21453038	21453038	CCATCTCATCCCTCGGTGTCCTCGACTCAGgtctggttcaaaACAAAGCTGGAAGACTCTTCCCTA	tctggttcaa
X	153587777	153587777	CCATCTCATCCCTCGGTGTCCTCGACTCAGgtctctataagCCAAACAGCCCAACAGTTTC	gtcctataag
X	153587777	153587777	CCATCTCATCCCTCGGTGTCCTCGACTCAGgtcagcctcGATGACCGGCTGTCTGTTT	gtcagcctcc
X	153587777	153587777	CCATCTCATCCCTCGGTGTCCTCGACTCAGtcaagctcgAAAGTGCCACCAACAA	tccaagctcg
19	44571260	44571260	CCATCTCATCCCTCGGTGTCCTCGACTCAGgttaaccagcgCTTGTAGCGCTTCCCACAGT	gtaccagcgc
19	44571260	44571260	CCATCTCATCCCTCGGTGTCCTCGACTCAGtctctattcggAGCTTCTTCCAAATCCTCA	tccattcggg
19	44571260	44571260	CCATCTCATCCCTCGGTGTCCTCGACTCAGgccagcgatctCTGTACCCCATAAATATGTACAA CT	gccagcgatt

TABLE 6C- continued

Chromosome	AlleleStart	AlleleEnd	Reverse	Barcode
X	153579431	153579431	CCATCTCATCCCTGCGTGTCTCCGACTCAGACTAGTAGTGGTAAGTGC	acctagactg
X	153579431	153579431	CCATCTCATCCCTGCGTGTCTCCGACTCAGACTGGTtcgctgcgcaaatcagtgacctcc	actggttcgc
X	153579431	153579431	CCATCTCATCCCTGCGTGTCTCCGACTCAGccaatattaggTCCCTTCCTGCGCACCTG	ccatattagg
12	46321441	46321441	CCATCTCATCCCTGCGTGTCTCCGACTCAGgctcgtcagcacaTGCTACTTTTGGGAATGAA G	gctcgtcagc
12	46321441	46321441	CCATCTCATCCCTGCGTGTCTCCGACTCAGgtaaatgacgCTTCTGAAACACCAATTGGAAA	cgtaaatgacg
12	46321441	46321441	CCATCTCATCCCTGCGTGTCTCCGACTCAGcggcgctgaTGTTAAGAGCCCAAGAGTTCA	ccggcgctga
X	153594210	153594210	CCATCTCATCCCTGCGTGTCTCCGACTCAGcgcgaagatacggggccctTACTCTTTGA	cgcaagata
X	153594210	153594210	CCATCTCATCCCTGCGTGTCTCCGACTCAGgaacccgcagactCGCAGCCCCCTACACTG	gaaccgcaga
X	153594210	153594210	CCATCTCATCCCTGCGTGTCTCCGACTCAGttggcagagataGACTGCCCTCTGCTGTG	ttggcagaga
16	3639306	3639306	CCATCTCATCCCTGCGTGTCTCCGACTCAGgcatctctgcAGTGACGATGAGCAGGAGGT	gcattctctgc
16	3639306	3639306	CCATCTCATCCCTGCGTGTCTCCGACTCAGttggaccgcagcCAATTCCCAATTGACCA	ttggaccgca
16	3639306	3639306	CCATCTCATCCCTGCGTGTCTCCGACTCAGggaagaacgtccCAAGCTTCTGAAACCAAGAC	gcagaacgtc
X	153599770	153599770	CCATCTCATCCCTGCGTGTCTCCGACTCAGaacttcgagcCTAGTGGGGGCATTCCAA	aacttcgagc
X	153599770	153599770	CCATCTCATCCCTGCGTGTCTCCGACTCAGgctcctagagCTTAGGGCGCGTTTCCT	gctcctagag
X	153599770	153599770	CCATCTCATCCCTGCGTGTCTCCGACTCAGtataetagtTTCAGCCTTTCCTCGCTCTA	tatctagctt
18	21453038	21453038	CCATCTCATCCCTGCGTGTCTCCGACTCAGgagatttggtTCCACATAACTCGCTTGCAG	gagtatgggc
18	21453038	21453038	CCATCTCATCCCTGCGTGTCTCCGACTCAGcctgagctcagAACTGTAGCCCAAGACATGC	cctgagctca
18	21453038	21453038	CCATCTCATCCCTGCGTGTCTCCGACTCAGcaggcgagtaAcaAAAGCTGGAAACTCTTCCCTA	caggcgagta
X	153587777	153587777	CCATCTCATCCCTGCGTGTCTCCGACTCAGcaggcagagCCAAGCCCAACAAGTTTC	gcaggcagag
X	153587777	153587777	CCATCTCATCCCTGCGTGTCTCCGACTCAGgctcgatatacGAATGACCGGCTGTCTGTTT	gcgtcgatac
X	153587777	153587777	CCATCTCATCCCTGCGTGTCTCCGACTCAGcagatgattatAAAGTGGCACCAACAACA	cgatgattat
19	44571260	44571260	CCATCTCATCCCTGCGTGTCTCCGACTCAGgacggctggcCTTGTACGCTTCCACACAGT	gacggctggc
19	44571260	44571260	CCATCTCATCCCTGCGTGTCTCCGACTCAGgagcctgagAGCTTCTTCCACAATCCTCA	ggagcctgag
19	44571260	44571260	CCATCTCATCCCTGCGTGTCTCCGACTCAGcctgactgctCTGTACCCCAATAAATGTACAACA CT	cctgactgct

Chromosome	AlleleStart	AlleleEnd	Reverse	Barcode
X	153579431	153579431	CCATCTCATCCCTGCGTGTC TCCGACTCAGAcggtcgaacgcGCCAGATGGGTAAAGTGC	aaggctgaag
X	153579431	153579431	CCATCTCATCCCTGCGTGTC TCCGACTCAGTaaccatagtTGCATAATCAGTGGCTCTCC	taaccatagc
X	153579431	153579431	CCATCTCATCCCTGCGTGTC TCCGACTCAGCtttgoccttcCTCCCTTCCTGCCACCTG	tettgccttc
12	46321441	46321441	CCATCTCATCCCTGCGTGTC TCCGACTCAGTtcttagattACATGTGTACTTTTGGGAATGAAG	tcttagatt
12	46321441	46321441	CCATCTCATCCCTGCGTGTC TCCGACTCAGTcatctcattCTTCTGAACA CCAAATTGGAAA	teatctcatt
12	46321441	46321441	CCATCTCATCCCTGCGTGTC TCCGACTCAGTctccgctsgTGTTAAGACC CCAGAGGTCA	tctccgctcg
X	153594210	153594210	CCATCTCATCCCTGCGTGTC TCCGACTCAGTgccatatgcggggcccCTACTCTTTGA	tgccatatgc
X	153594210	153594210	CCATCTCATCCCTGCGTGTC TCCGACTCAGTaaggcctctCTCGAGCCCCCTACACTG	taaggcctct
X	153594210	153594210	CCATCTCATCCCTGCGTGTC TCCGACTCAGAgagtaggcogTGACTGCCCTCTGCTGTG	gagt aggcog
16	3639306	3639306	CCATCTCATCCCTGCGTGTC TCCGACTCAGGcaaat aagctAGTGACGATGAGCAGGAGGT	gcaataagct
16	3639306	3639306	CCATCTCATCCCTGCGTGTC TCCGACTCAGBgpggt tgcaggCCAATTC CCAATTGACCA	ggcgttgcaa
16	3639306	3639306	CCATCTCATCCCTGCGTGTC TCCGACTCAGccaagaagcgCCAAGTTCTTGAACACAGAC	cgaagaagcg
X	153599770	153599770	CCATCTCATCCCTGCGTGTC TCCGACTCAGbggttaacctcgCTAGTGGGGGCAATTC CAA	ggttacctcg
X	153599770	153599770	CCATCTCATCCCTGCGTGTC TCCGACTCAGctccgccttaCTAGGGGCGCGTTTCTCT	ctccgcctta
X	153599770	153599770	CCATCTCATCCCTGCGTGTC TCCGACTCAGctccagagatTCAGCCTTTTCCTCGCTCTA	ctccagagat
18	21453038	21453038	CCATCTCATCCCTGCGTGTC TCCGACTCAGgtcgaggt agTCCACATAACTCGCTTG CAG	gtcgaggt ag
18	21453038	21453038	CCATCTCATCCCTGCGTGTC TCCGACTCAGTatggacctggAACTGTAGCC CCAAGACACTGC	tatggacctg
18	21453038	21453038	CCATCTCATCCCTGCGTGTC TCCGACTCAGTacctgctagACAAAAGCTGG AAACTCTTCCCTA	tacctgctag
X	153587777	153587777	CCATCTCATCCCTGCGTGTC TCCGACTCAGcgcgcgcgcgcgCAACAAGCCC AACAAGTTC	cgcgcgcgcg
X	153587777	153587777	CCATCTCATCCCTGCGTGTC TCCGACTCAGgttgaacgttGAATGACCCGGCTGTCTGTTT	gttgaacgtt
X	153587777	153587777	CCATCTCATCCCTGCGTGTC TCCGACTCAGTgcacaacgcAAAAGTGGCAC CCAACAA	tgccaacgca
19	44571260	44571260	CCATCTCATCCCTGCGTGTC TCCGACTCAGggattgacctCTTGTAGCGCTTCC CACAGT	ggattgacct
19	44571260	44571260	CCATCTCATCCCTGCGTGTC TCCGACTCAGBgacggattcAGCTTCTTTCCACAATCTCA	ggacggattc
19	44571260	44571260	CCATCTCATCCCTGCGTGTC TCCGACTCAGTcctcsgtsgCTGTACCC CCAATAATGTACAACA CT	tctccgctcg
19	153579431	153579431	CCATCTCATCCCTGCGTGTC TCCGACTCAGagtttcaggtCCGCCAGATGGGTAAAGTGC	agttcatggt

TABLE 6C- continued

Chromosome	AlleleStart	AlleleEnd	Reverse	Barcode
X	153579431	153579431	CCATCTCATCCCTGCGTGTCTCCGACTCAGTatccattccTGCAAAATCAGTGGCTCTCC	tatccattcc
X	153579431	153579431	CCATCTCATCCCTGCGTGTCTCCGACTCAGGagagcgcgCTCCCTTCTCTGCCACCTG	ggagagcgcg
12	46321441	46321441	CCATCTCATCCCTGCGTGTCTCCGACTCAGggaccttggACATGTGATACTTTTGGGAATGAA G	cggaccttgg
12	46321441	46321441	CCATCTCATCCCTGCGTGTCTCCGACTCAGggcaatctccctTCTGAAACACCAATTGGAAA	ggcaatctcc
12	46321441	46321441	CCATCTCATCCCTGCGTGTCTCCGACTCAGaggtattgatTTGTTAAGAGCCAGAGTTTCA	aggattgatt
X	153594210	153594210	CCATCTCATCCCTGCGTGTCTCCGACTCAGgcegttgccctGGGGCCCCCTACTCTTTGA	gccgttgccct
X	153594210	153594210	CCATCTCATCCCTGCGTGTCTCCGACTCAGAagtaacgtcgtCTCGCAGCCCCCTACACTG	aagtaacgtcg
X	153594210	153594210	CCATCTCATCCCTGCGTGTCTCCGACTCAGtggcttaaggtGACTGCCCTCTCTCTGTG	tggtttaagg
16	3639306	3639306	CCATCTCATCCCTGCGTGTCTCCGACTCAGctcttceagaAGTGACGATGAGCAGGAGGT	ctcttccaga
16	3639306	3639306	CCATCTCATCCCTGCGTGTCTCCGACTCAGegttcttcaagCCAAATTCCCAATFGACCA	cgttcttcaa
16	3639306	3639306	CCATCTCATCCCTGCGTGTCTCCGACTCAGcaagcgctgcCCAAGCTTCTGAAACCAGAC	caacggctgc
X	153599770	153599770	CCANCTCATCCCTGCGTGTCTCCGACTCAGgcaagtaaacctAGTGGGGCATTCCAA	gcaagtaaac
X	153599770	153599770	CCATCTCATCCCTGCGTGTCTCCGACTCAGgttcataagtcCTCTAGGGCGCGTTTCTT	gttcataagtc
X	153599770	153599770	CCANCTCATCCCTGCGTGTCTCCGACTCAGacggcgagccctCAGCCCTTCTCTCGCTCTTA	acggcgagcc
18	21453038	21453038	CCATCTCATCCCTGCGTGTCTCCGACTCAGgtatgggtcggTCCACATAACTCGCTTGCAG	gtatggtcgg
18	21453038	21453038	CCATCTCATCCCTGCGTGTCTCCGACTCAGtccggttatccgAACTGTAGCCCAAGACATGC	tccggttatcc
18	21453038	21453038	CCATCTCATCCCTGCGTGTCTCCGACTCAGggggtcgataaCAAAAGCTGGAAGACTTTCCCTA	gcggtcgata
X	153587777	153587777	CCATCTCATCCCTGCGTGTCTCCGACTCAGtccctcagtatCCAACAAGCCCAACAAGTTTC	tccctcagtat
X	153587777	153587777	CCATCTCATCCCTGCGTGTCTCCGACTCAGaccgttccctgGAATGACCGGCTGTCTGTTT	accgttccctg
X	153587777	153587777	CCATCTCATCCCTGCGTGTCTCCGACTCAGgctcgtcttAAAGTGGCACCAACAACA	gccctgctctt
19	44571260	44571260	CCATCTCATCCCTGCGTGTCTCCGACTCAGagcgtaaccaCTTGTAGCGGTTTCCCACAGT	agcgtaacca
19	44571260	44571260	CCATCTCATCCCTGCGTGTCTCCGACTCAGtggctgatgAGCTTCTTTCACAAATCCTCA	ttgctgatg
19	44571260	44571260	CCANCTCATCCCTGCGTGTCTCCGACTCAGtattgatctCTGTACCCCAATAAATGTACAACA CT	ttattgatct
X	153579431	153579431	CCATCTCATCCCTGCGTGTCTCCGACTCAGtaccgcteggaCGCCAGATGGGTAAATGTC	taccgctegga
X	153579431	153579431	CCATCTCATCCCTGCGTGTCTCCGACTCAGcaaatccaaggtGCAAAATCAGTGGCTCTCC	caatccaagg

TABLE 6C- continued

Chromosome	AlleleStart	AlleleEnd	Reverse	Barcode
X	153579431	153579431	CCATCTCATCCCTGCGTGTCTCCGACTCAGTcgtagetatCTCCCTTCCTGCGCACCTG	tcgtagetat
12	46321441	46321441	CCATCTCATCCCTGCGTGTCTCCGACTCAGcgctcatcgACAATGTATACTTTTGGGAATGAA G	cgctcatcgc
12	46321441	46321441	CCATCTCATCCCTGCGTGTCTCCGACTCAGTccggtteattCTTCTGAAACACCAAAATTGGAAA	tccggttcatt
12	46321441	46321441	CCATCTCATCCCTGCGTGTCTCCGACTCAGggccaggctTGTTAAGAGCCAGAGGTTCA	cgggccaggct
X	153594210	153594210	CCATCTCATCCCTGCGTGTCTCCGACTCAGaaacctatctGGGGCCCCCTACTCTTTGA	caacctatct
X	153594210	153594210	CCATCTCATCCCTGCGTGTCTCCGACTCAGegtaaatetcaCTCGACGCCCTACACTG	cgtaaatetca
X	153594210	153594210	CCATCTCATCCCTGCGTGTCTCCGACTCAGatatcgcgactGACTGCCCCCTCTGCTGTG	atatcgcgac
16	3639306	3639306	CCATCTCATCCCTGCGTGTCTCCGACTCAGTcaatatctgAGTGACGATGAGCAGGAGGT	tcaatatctg
16	3639306	3639306	CCATCTCATCCCTGCGTGTCTCCGACTCAGtagagtagatagCCAAATTCCTCCATTGACCA	atagagtagata
16	3639306	3639306	CCATCTCATCCCTGCGTGTCTCCGACTCAGgcaactagttCCAAGCTTCCTGAACCAAGAC	gcaactagtt
X	153599770	153599770	CCATCTCATCCCTGCGTGTCTCCGACTCAGatctogaaatcCTAGTGGGGGCATTCCAA	atctcgaatc
X	153599770	153599770	CCATCTCATCCCTGCGTGTCTCCGACTCAGccaggagcgactCTAGGGCGCGTTTCCT	ccaggagcgga
X	153599770	153599770	CCATCTCATCCCTGCGTGTCTCCGACTCAGatctccatcgTCAGCCCTTTCTCGCTCTTA	atctccatcg
18	21453038	21453038	CCATCTCATCCCTGCGTGTCTCCGACTCAGttgacgagctTCCACATAACTCGCTTGCAG	ttgacgagct
18	21453038	21453038	CCATCTCATCCCTGCGTGTCTCCGACTCAGtactattaccGAAGTGTAGCCCGAGACATGTC	tactattacc
18	21453038	21453038	CCATCTCATCCCTGCGTGTCTCCGACTCAGcgtcctgggacACAAAAGCTGGAAACTCTTCCCTA	cgctcctggac
X	153587777	153587777	CCATCTCATCCCTGCGTGTCTCCGACTCAGetcgggcgcttCCAACAGCCCAACAAGTTC	ctcggcgctt
X	153587777	153587777	CCATCTCATCCCTGCGTGTCTCCGACTCAGgat acgtaaGATGACCGGCTGTCTGTTTT	gatacgtaag
X	153587777	153587777	CCATCTCATCCCTGCGTGTCTCCGACTCAGetcgsgattaaAAAAGTGGCACCAACAACA	ctcggattaa
19	44571260	44571260	CCATCTCATCCCTGCGTGTCTCCGACTCAGttggattcgtCTTGTAGCGGTTCCACACAGT	ttggattcgt
19	44571260	44571260	CCATCTCATCCCTGCGTGTCTCCGACTCAGcgtccgctaaAGCTTCTTTTCCACAATCCTCA	ccgtccgcta
19	44571260	44571260	CCATCTCATCCCTGCGTGTCTCCGACTCAGgggattgcaactGTATACCCCAATAAATGTACAAAC ACT	gcgattgcaa
X	153579431	153579431	CCATCTCATCCCTGCGTGTCTCCGACTCAGccatgcataaCGCCAGATGGGTAAAGTGC	ccatgcataa
X	153579431	153579431	CCATCTCATCCCTGCGTGTCTCCGACTCAGtaattgcaatTGCAAAATCAGTGGCTCTCC	taattgcaat

TABLE 6C- continued

Chromosome	AlleleStart	AlleleEnd	Reverse	Barcode
X	153579431	153579431	CCATCTCATCCCTCGGTGTCCTCGACTCAGacgactccaaCTCCTTCCTGCGCACCTG	acgactccaa
12	46321441	46321441	CCATCTCATCCCTCGGTGTCCTCGACTCAGatcatgagaACATGTGATACTTTTGGGAATGAA G	atcatgacga
12	46321441	46321441	CCATCTCATCCCTCGGTGTCCTCGACTCAGaactcctaactCTTCTGAACACCAAAATTGGAAA	aactcctaact
12	46321441	46321441	CCATCTCATCCCTCGGTGTCCTCGACTCAGggatattogtTGTTAAGAGCCCCAGAGGTTCA	ggatattogt
X	153594210	153594210	CCATCTCATCCCTCGGTGTCCTCGACTCAGTcggatgactgggccccCTACTCTTTGA	tggatgact
X	153594210	153594210	CCATCTCATCCCTCGGTGTCCTCGACTCAGgacgogegagCTCGACGCCCTACACTG	gacgogegag
X	153594210	153594210	CCATCTCATCCCTCGGTGTCCTCGACTCAGgctagacctTGACTGCCCTCTGCTGTG	gcttagacct
16	3639306	3639306	CCATCTCATCCCTCGGTGTCCTCGACTCAGgacaggaagAGTAGACATGAGCAGAGGT	gaccagggaga
16	3639306	3639306	CCATCTCATCCCTCGGTGTCCTCGACTCAGgctetggcgtGCCAATTCCCATTGACCA	gctctggcgt
16	3639306	3639306	CCATCTCATCCCTCGGTGTCCTCGACTCAGTgggtccggaaCCAAAGCTTCGTGAACCAGAC	tggtccggaa
X	153599770	153599770	CCATCTCATCCCTCGGTGTCCTCGACTCAGctctgctctCTAGTGGGGCATTCCAA	ctctgctct
X	153599770	153599770	CCATCTCATCCCTCGGTGTCCTCGACTCAGcagaagcagCTTAGGGCGCGTTTCCT	ccagaagcag
X	153599770	153599770	CCATCTCATCCCTCGGTGTCCTCGACTCAGggaaggttgctCAGCCTTTTCCTCGCTCTA	ggaaggttgc
18	21453038	21453038	CCATCTCATCCCTCGGTGTCCTCGACTCAGtaacgggtacgTCCACATPACTCGCTTGCAAG	taacgggtacg
18	21453038	21453038	CCATCTCATCCCTCGGTGTCCTCGACTCAGctcgtcteatgGAAGTGTAGCCCCAGACACTGC	ctcgtcteatg
18	21453038	21453038	CCATCTCATCCCTCGGTGTCCTCGACTCAGactccaaaggcACAAAGCTGGAAACTCTTCCCTA	actccaaaggc
X	153587777	153587777	CCATCTCATCCCTCGGTGTCCTCGACTCAGgagctgetatCCAAACAGCCCCAACAGTTTC	gagctgetat
X	153587777	153587777	CCATCTCATCCCTCGGTGTCCTCGACTCAGcgttgaggccGAATGACCGGCTGCTGTGTTT	cgttgaggcc
X	153587777	153587777	CCATCTCATCCCTCGGTGTCCTCGACTCAGTctctggatccAAAGTGCACACCAACAA	ttctggatcc
19	44571260	44571260	CCATCTCATCCCTCGGTGTCCTCGACTCAGcgggattccactTGTAGCGCTTCCCACAGT	ccggattccca
19	44571260	44571260	CCATCTCATCCCTCGGTGTCCTCGACTCAGTccatcgcttAGCTTCTTTCCCAATTCCTCA	tccatcgctt
19	44571260	44571260	CCATCTCATCCCTCGGTGTCCTCGACTCAGTtaactctcaCTGTACCCCATAAATATGTACAACA CT	ttacttctca

Sensitivity and Reproducibility of Assay

[0164] The AAF of somatic mutations can vary dramatically across tissues, where they can be nearly undetectable in tissues such as blood, but higher frequency in tissues like the brain. Given that most genetic testing is performed on blood or cell free DNA samples with anticipated low AAFs, the ability of the presently described methods to accurately detect AAFs at extremely low levels, which are often difficult or impossible to accurately assess by other methods.

[0165] The sensitivity of triple-primer PCR sequencing was assessed through serial dilution of a genomic control DNA sample containing the same 5 known germline mutations described above (Tables 6A-6C) with a control DNA lacking these mutations, thereby generating AAFs ranging from 50% down to 0.01%. The dilutions were amplified with primers for each mutation and sequenced on the Ion Torrent S5 with sequencing reads of 400 bp in length. All reads were processed using custom analytical scripts (described in methods), allowing the comparison of assessed and expected allelic fractions.

[0166] The presently described method accurately measures AAFs as low as 0.01% when using a 50 ng of genomic DNA, although for significant detection above the amplicon-specific error rates, AAFs were typically required to be above 0.05% (FIGS. 18A, 18B). Surprisingly, 6 of 6 mutations were successfully identified at AAFs of 0.05%, and all were identified by at least one of the primers in the sets at AAFs as low as 0.01%. Therefore, the presently described approach is able to achieve a 100% sensitivity for detection of alleles down to 0.01% AAF (FIGS. 18A, 18B). The largest factors observed in accurately measuring the AAFs at extremely low levels of below 0.05% was providing sufficient input DNA and achieving enough sequencing depth to distinguish errors from true calls. In this case, a depth of more than 50,000× is recommended for the best sensitivity. While each independent primer set can produce slightly different AAFs due to both inherent primer characteristics and variability amongst reactions, averaging the primers provides an extremely accurate assessment of the true AAF. Even more, the accuracy of the estimate is better assessed through the comparison of the confidence intervals from the AAFs of the mutation and the background error rates. For example, it was found that the measurement of a 2048-fold dilution (estimated AAF~0.012%) sample resulted in an AAF of 0.0136%±0.012% while the background error rate was significantly lower than the measured AAF at 0.0015%±0.009%.

[0167] The measured AAFs (average across triple primer sets) were linearly correlated with the expected AAFs down to 0.01% ($R^2>0.999$), though as expected, individual AAFs do vary amongst individual primers ($R^2>0.98$). Therefore, while individual primer sets are prone to biases in AAFs, the utility of multiple primer provides a robust and accurate measurement.

[0168] DNA is often limited, particularly in clinical contexts, but is also known as an important factor for sensitivity for somatic alleles due to the presence of fewer DNA fragments containing the targeted allele. Therefore, the sensitivity of using 50 ng was compared to using a reduced concentration of 25 ng (~3800 cells) (PMID: 30813969). With 3800 cells, the accurate detection of the lowest dilution of 0.01% AAF is unlikely as it would likely only be represented by a single fragment. Surprisingly, AAFs down to 0.05% remained detectable with 25 ng DNA (FIGS. 18C,

18D), though with less precision, which indicates that increasing the input DNA to 50 ng or more would improve accuracy when validating alleles below 0.1% AAF.

[0169] Furthermore, the impact of total sequencing depth on the accuracy was assessed to identify the minimum depth needed for accurate determination of AAFs. Sequencing data for each amplicon were randomly sampled to create artificial datasets containing a wide range of depths ranging from 10,000 to 150,000× coverage. Increasing read depths above 10,000× did not have a substantial impact on the background error rates within the amplicons. Even more, a minimum depth of 10,000× was able to accurately measure AAFs down to 0.1% with no improvement with elevated coverage. However, accurate measurement of AAFs below 0.1% required depths of 25,000× to ensure significance over the background errors. Overall, a strong correlation was found of AAFs measured across a wide range of read depths, indicating that detection of AAFs of 0.01% is possible at depths greater than above 25,000×.

[0170] The assessment of error rates and the potential for false positive allele calls was extended by performing similar sequencing on DNA samples lacking mutations. As expected, these alleles were not detectable, with only the typical background error rate being detected, which is often not the same allele as the mutation, supporting the specificity of this method.

Precise Assessment of Broad Range of AAFs in Multiple Tissues

[0171] As some tissues are more difficult to work with, the ability was assessed of the method to accurately detect known mosaic alleles that were previously identified in blood and brain tissue by a range of methods including WGS, WES, and targeted Illumina sequencing. Even more, given the importance of validating indels and the elevated indels error rates on Ion Torrent data, >50 somatic indels were tested using the method of the present invention with a direct comparison of the sites between the DNA sample containing the mutation and a control sample. It was demonstrated that AAFs of SNVs ($R=0.93$, (FIG. 17A) and indels ($R=0.89$, across insertions and deletions (FIGS. 19A, 19B)) detected between the methods were highly correlated regardless of the tissue or original sequencing platform. Surprisingly, very accurate assessments of indels with very little increase in error rates were obtained. However, the ability to validate extremely low AAF indels occurring within homopolymers remained challenging when using Ion Torrent. In some instances, AAFs were observed that were dissimilar to the original detection method. In these instances, the discrepancy was driven by low coverage in the original sequencing platform, resulting in an incorrect estimate of AAFs. Additionally, in some cases, a single primer provided an outlier AAF, which deviated from the other primers and original method of identification. In these cases, other primers revealed a germline mutation impacting the primer binding, resulting in allelic dropout. Such instances of allelic dropout are mitigated through the primer design process, but as often is the case, not all alleles are known, particularly in targeted sequencing and exome studies. The chances of allelic dropout highlight the importance of using multiple primers when studying mosaic and germline alleles.

Robust Validation for Low AAF Insertions/Deletions

[0172] The known increased error rates for indel in Ion Torrent data and the inability to utilize PCR duplicate information may limit the ability to quantitate some ultra-rare alleles (<0.05% AAF) and indels. Even more, the Pollux software is known to overcorrect for indels and has difficulty distinguishing rare indels from artifacts. Despite these limitations, it was assessed how the method performs on a wide range of indels occurring at AAFs from 1% to 30% and 1 to 21 base pairs in length, including 40 insertions and 60 deletions previously identified using 200× whole genome sequencing. Even more importantly, these mutations were not identified in control DNA, where at these sites very low error rates for indels ($0.010\% \pm 0.05\%$) were found, supporting that even the single base indels are not being introduced by PCR or the Ion Torrent. These data indicate a sensitivity to accurately quantitate AAFs of indels down to 0.05% in many instances. Despite that many of these mutations were detected using only a few reads in the WGS data, a strong correlation was found between the predicted AAFs in the WGS and the measured values by the method described in this example (FIGS. 19A, 19B; $R^2=0.75$ deletions and $R^2=0.94$ for insertions), indicating that this method is also sensitive to detect very low AAF indels, which are often difficult to validate.

[0173] To further improve the sensitivity for low AAFs, a modified version of the protocol was performed (FIG. 5A) in which an initial low cycle PCR was performed containing biotinylated dCTP (~25% of a cytosines) and using unique molecular indexes (UMIs) to uniquely tag all PCR products in the first 10 cycles. After purification using either streptavidin capture or enzymatic digest (see methods), all reactions were further amplified by a common primer that maintained the UMI signature, effectively tagging all PCR duplicates from the second round of PCR. An optional step after purification comprises analyzing the sample for acceptable quality control, which, for example, can be done using a Bioanalyzer or TapeStation (FIG. 5B).

[0174] The incorporation of biotin into the PCR product did not impact the overall measured AAFs, but slightly reduced the error rate ($0.0023\% \pm 0.0011\%$ AAF), possibly due to the ability to perform better purification and the use of a common primer for the majority of the amplifications. These indicate that a 2-step UMI approach for the method is valuable in situations requiring reduced error rates for ultra-low AAFs or where PCR duplicates may be of particular concern.

Application of Method for Novel Variant Discovery Using Illumina Sequencing

[0175] The increased sensitivity of the the presently described approach can be further applied for the detection of novel ultra-low AAFs variants with Illumina-based sequencing. Overlapping primers were developed so that all regions of the PRNP gene was covered by at least 3 independent amplicons, each containing Illumina sequencing adapters and UMIs. Using the 2-step PCR approach, sequencing libraries were prepared for a dilution series of a known mutation (5%, 0.5%, and 0.05% AAFs) and additional samples were screened for novel alleles. While any given amplicon can have some errors, as outlined above and previously documented in amplicon-based sequencing studies, it was contemplated whether the method could reduce

such effects to identify high-confidence mutations. By requiring consistent AAFs across multiple unique primer sets, the AAFs of mutations were accurately measured down to at least 0.05% (FIG. 19C). Even more, when applied to a large set of tissues derived DNA samples for detections of novel mutations in a given gene, mutations down to 0.05% AAFs were accurately detected with no additional false positive occurrences (FIGS. 19C and 19D), indicating a possible option for improved accurate measurement of AAFs of novel alleles in targeted sequencing platforms.

[0176] The following materials and methods were used in carrying out this example.

Primer Design

[0177] At least three unique sets of primers were designed for each mutation by extracting the flanking sequence around each mutation so that the mutation is located at different positions within each of the three sequences. Next, common alleles are masked, along with the targeted mutation and flanking 5bps on each site using the bedtools maskfasta tool. The masked multi-fasta file containing all sequences for targeted alleles are input into BatchPrimer webtool to design primers for each sequence. Primers are designed to an average TM of 60° C., with a minimum of 59° C. and maximum of 62° C. The amplicon length is dependent on the specific mutation and DNA sources. For example, difficult to map regions may have longer products while degraded DNA samples may require shorter amplicons. In general, to ensure that all primers are likely unique and of similar amplicon length, amplicons have a target length of 225-300 bp in length. The primer sequences are checked by BLAT and in-silico PCR to ensure both their unique amplicon in the genome and that the primer binding sites do not overlap between any set of primers. The final set of primers are then uniquely barcoded using 10 nt barcodes and if desired, an additional 10 nt UMI is added. Finally, Ion Torrent specific adapter sequences are appended to the forward and reverse primers, allowing for their direct sequencing.

Library Preparation

[0178] For the standard, single step PCR sequencing method described above, PCR was performed using 20 cycles on a 25 µl reaction mix containing either 25 or 50 ng of input DNA sample, Phusion Hot-Start polymerase, dNTPs, HC-Buffer, and the primers. For initial testing, 30 cycles of enrichment were used to ensure only a single amplicon is produced. The high-sensitivity method modifies this process by reduction of the PCR cycling to 5 and the incorporation of 0.1 µL of 0.4 mM biotin-14-dCTP into the reaction mix. Biotinylated PCR amplicons are captured by adding 5 µl of washed Streptavidin Myone beads resuspended in 25 µl of 2× binding and washing buffer. The mixture is incubated at room temperature with gentle mixing for 15 minutes and placed on a 96-well magnetic plate. The liquid was removed and the beads were washed one time with 1× binding and washing buffer. Then beads are then resuspended in 25 µl PCR reaction mixture containing custom primers which preserve the original UMI sequences, Phusion Hot-Start polymerase, dNTPs, and HC-Buffer. The biotin labeled product was amplified with an additional 20 cycles of enrichment before the beads were removed.

Enriched products were pools at equal volumes and purified using the MagJet purification kit.

QC and Variant Calling

[0179] Purified library pools are analyzed for enrichment efficiency and the complete removal of primers through by either the Agilent Bioanalyzer Hi-sensitivity chip or the TapeStation. The concentration was determined using PicoGreen. Pools were diluted to a final concentration of 100 pM prior to sequencing on the 430 chip for the Ion Torrent S5.

[0180] Raw unmapped bam files were obtained for each run and were processed using our custom analyses pipeline. First, all BAMs are converted to a fastq file using bedtools bamtofastq tool. Then, quality and adapter trimming was performed using cutadapt tool. Next, samples lacking UMIs, are demultiplexed using fastx_barcode_splitter, resulting in separate fastq files for each primer set. The barcode sequences are removed from the sequences using cutadapt. If the allele being tested in an SNV, indel correction is performed using Pollux. Finally, all samples are aligned to the reference genome using BWA-mem.

[0181] Variants are then called across the length of each amplicon through the use of samtools mPileup with the

settings: q=20, Q=20. The resulting vcfs are parsed into a file containing the flanking 50 nt positions on each side of the variant and a separate file for the allele of interest. The average allele frequency across the flanking regions are then compared to the average AAF of the mutation across the 3 unique primers.

OTHER EMBODIMENTS

[0182] From the foregoing description, it will be apparent that variations and modifications may be made to the invention described herein to adopt it to various usages and conditions. Such embodiments are also within the scope of the following claims.

[0183] The recitation of a listing of elements in any definition of a variable herein includes definitions of that variable as any single element or combination (or subcombination) of listed elements. The recitation of an embodiment herein includes that embodiment as any single embodiment or in combination with any other embodiments or portions thereof.

[0184] All patents and publications mentioned in this specification are herein incorporated by reference to the same extent as if each independent patent and publication was specifically and individually indicated to be incorporated by reference.

SEQUENCE LISTING

<160> NUMBER OF SEQ ID NOS: 339

<210> SEQ ID NO 1

<211> LENGTH: 18

<212> TYPE: DNA

<213> ORGANISM: Artificial Sequence

<220> FEATURE:

<223> OTHER INFORMATION: Description of Artificial Sequence: Synthetic primer

<400> SEQUENCE: 1

cagggcctca ccttggtc

18

<210> SEQ ID NO 2

<211> LENGTH: 22

<212> TYPE: DNA

<213> ORGANISM: Artificial Sequence

<220> FEATURE:

<223> OTHER INFORMATION: Description of Artificial Sequence: Synthetic primer

<400> SEQUENCE: 2

ctgtgacata gcactcctcc ag

22

<210> SEQ ID NO 3

<211> LENGTH: 18

<212> TYPE: DNA

<213> ORGANISM: Artificial Sequence

<220> FEATURE:

<223> OTHER INFORMATION: Description of Artificial Sequence: Synthetic primer

<400> SEQUENCE: 3

aggctggctg gttgacct

18

<210> SEQ ID NO 4

<211> LENGTH: 26

-continued

<212> TYPE: DNA
<213> ORGANISM: Artificial Sequence
<220> FEATURE:
<223> OTHER INFORMATION: Description of Artificial Sequence: Synthetic primer

<400> SEQUENCE: 4

aatcacactc cataggtatc atttca 26

<210> SEQ ID NO 5
<211> LENGTH: 25
<212> TYPE: DNA
<213> ORGANISM: Artificial Sequence
<220> FEATURE:
<223> OTHER INFORMATION: Description of Artificial Sequence: Synthetic primer

<400> SEQUENCE: 5

ttcattcatt tgtttaagat cagca 25

<210> SEQ ID NO 6
<211> LENGTH: 28
<212> TYPE: DNA
<213> ORGANISM: Artificial Sequence
<220> FEATURE:
<223> OTHER INFORMATION: Description of Artificial Sequence: Synthetic primer

<400> SEQUENCE: 6

ttaacggacg cgccagatgg gtaagtgc 28

<210> SEQ ID NO 7
<211> LENGTH: 29
<212> TYPE: DNA
<213> ORGANISM: Artificial Sequence
<220> FEATURE:
<223> OTHER INFORMATION: Description of Artificial Sequence: Synthetic primer

<400> SEQUENCE: 7

tccggcttac tgcaaatcag tggctctcc 29

<210> SEQ ID NO 8
<211> LENGTH: 28
<212> TYPE: DNA
<213> ORGANISM: Artificial Sequence
<220> FEATURE:
<223> OTHER INFORMATION: Description of Artificial Sequence: Synthetic primer

<400> SEQUENCE: 8

tctcattcag ctcccttctc gccacctg 28

<210> SEQ ID NO 9
<211> LENGTH: 35
<212> TYPE: DNA
<213> ORGANISM: Artificial Sequence
<220> FEATURE:
<223> OTHER INFORMATION: Description of Artificial Sequence: Synthetic primer

<400> SEQUENCE: 9

gcggtcatac acatgtgata cttttgggaa tgaag 35

-continued

<210> SEQ ID NO 10
<211> LENGTH: 32
<212> TYPE: DNA
<213> ORGANISM: Artificial Sequence
<220> FEATURE:
<223> OTHER INFORMATION: Description of Artificial Sequence: Synthetic primer

<400> SEQUENCE: 10
taggacgttc cttctgaaca ccaaattgga aa 32

<210> SEQ ID NO 11
<211> LENGTH: 15
<212> TYPE: DNA
<213> ORGANISM: Artificial Sequence
<220> FEATURE:
<223> OTHER INFORMATION: Description of Artificial Sequence: Synthetic oligonucleotide

<400> SEQUENCE: 11
ttaacggacg cgcca 15

<210> SEQ ID NO 12
<211> LENGTH: 15
<212> TYPE: DNA
<213> ORGANISM: Artificial Sequence
<220> FEATURE:
<223> OTHER INFORMATION: Description of Artificial Sequence: Synthetic oligonucleotide

<400> SEQUENCE: 12
tccggcttac tgcaa 15

<210> SEQ ID NO 13
<211> LENGTH: 15
<212> TYPE: DNA
<213> ORGANISM: Artificial Sequence
<220> FEATURE:
<223> OTHER INFORMATION: Description of Artificial Sequence: Synthetic oligonucleotide

<400> SEQUENCE: 13
tctcattcag ctccc 15

<210> SEQ ID NO 14
<211> LENGTH: 15
<212> TYPE: DNA
<213> ORGANISM: Artificial Sequence
<220> FEATURE:
<223> OTHER INFORMATION: Description of Artificial Sequence: Synthetic oligonucleotide

<400> SEQUENCE: 14
gcggtcatac acatg 15

<210> SEQ ID NO 15
<211> LENGTH: 15
<212> TYPE: DNA
<213> ORGANISM: Artificial Sequence
<220> FEATURE:
<223> OTHER INFORMATION: Description of Artificial Sequence: Synthetic oligonucleotide

<400> SEQUENCE: 15
taggacgttc cttct 15

-continued

<210> SEQ ID NO 16
<211> LENGTH: 15
<212> TYPE: DNA
<213> ORGANISM: Artificial Sequence
<220> FEATURE:
<223> OTHER INFORMATION: Description of Artificial Sequence: Synthetic
oligonucleotide

<400> SEQUENCE: 16

caaggtgagg ccctg 15

<210> SEQ ID NO 17
<211> LENGTH: 15
<212> TYPE: DNA
<213> ORGANISM: Artificial Sequence
<220> FEATURE:
<223> OTHER INFORMATION: Description of Artificial Sequence: Synthetic
oligonucleotide

<400> SEQUENCE: 17

agtgctatgt cacag 15

<210> SEQ ID NO 18
<211> LENGTH: 15
<212> TYPE: DNA
<213> ORGANISM: Artificial Sequence
<220> FEATURE:
<223> OTHER INFORMATION: Description of Artificial Sequence: Synthetic
oligonucleotide

<400> SEQUENCE: 18

tcaaccagcc agcct 15

<210> SEQ ID NO 19
<211> LENGTH: 15
<212> TYPE: DNA
<213> ORGANISM: Artificial Sequence
<220> FEATURE:
<223> OTHER INFORMATION: Description of Artificial Sequence: Synthetic
oligonucleotide

<400> SEQUENCE: 19

ctatggagtg tgatt 15

<210> SEQ ID NO 20
<211> LENGTH: 15
<212> TYPE: DNA
<213> ORGANISM: Artificial Sequence
<220> FEATURE:
<223> OTHER INFORMATION: Description of Artificial Sequence: Synthetic
oligonucleotide

<400> SEQUENCE: 20

aaacaaatga atgaa 15

<210> SEQ ID NO 21
<211> LENGTH: 41
<212> TYPE: DNA
<213> ORGANISM: Artificial Sequence
<220> FEATURE:
<223> OTHER INFORMATION: Description of Artificial Sequence: Synthetic
primer

<400> SEQUENCE: 21

-continued

cctctctatg ggcagtcggt gatcagggcc tcaccttggt c 41

<210> SEQ ID NO 22
<211> LENGTH: 45
<212> TYPE: DNA
<213> ORGANISM: Artificial Sequence
<220> FEATURE:
<223> OTHER INFORMATION: Description of Artificial Sequence: Synthetic
primer

<400> SEQUENCE: 22

cctctctatg ggcagtcggt gatctgtgac atagcactcc tccag 45

<210> SEQ ID NO 23
<211> LENGTH: 41
<212> TYPE: DNA
<213> ORGANISM: Artificial Sequence
<220> FEATURE:
<223> OTHER INFORMATION: Description of Artificial Sequence: Synthetic
primer

<400> SEQUENCE: 23

cctctctatg ggcagtcggt gataggtctg ctggttgacc t 41

<210> SEQ ID NO 24
<211> LENGTH: 49
<212> TYPE: DNA
<213> ORGANISM: Artificial Sequence
<220> FEATURE:
<223> OTHER INFORMATION: Description of Artificial Sequence: Synthetic
primer

<400> SEQUENCE: 24

cctctctatg ggcagtcggt gataatcaca ctccataggt atcatttca 49

<210> SEQ ID NO 25
<211> LENGTH: 48
<212> TYPE: DNA
<213> ORGANISM: Artificial Sequence
<220> FEATURE:
<223> OTHER INFORMATION: Description of Artificial Sequence: Synthetic
primer

<400> SEQUENCE: 25

cctctctatg ggcagtcggt gatttcattc atttgtttaa gatcagca 48

<210> SEQ ID NO 26
<211> LENGTH: 46
<212> TYPE: DNA
<213> ORGANISM: Artificial Sequence
<220> FEATURE:
<223> OTHER INFORMATION: Description of Artificial Sequence: Synthetic
primer

<400> SEQUENCE: 26

cctctctatg ggcagtcggt gattcaatgt gtgttttagg caactc 46

<210> SEQ ID NO 27
<211> LENGTH: 41
<212> TYPE: DNA
<213> ORGANISM: Artificial Sequence
<220> FEATURE:
<223> OTHER INFORMATION: Description of Artificial Sequence: Synthetic
primer

-continued

<400> SEQUENCE: 27

cctctctatg ggcagtcggt gataggggga catgcaagac a 41

<210> SEQ ID NO 28

<211> LENGTH: 41

<212> TYPE: DNA

<213> ORGANISM: Artificial Sequence

<220> FEATURE:

<223> OTHER INFORMATION: Description of Artificial Sequence: Synthetic primer

<400> SEQUENCE: 28

cctctctatg ggcagtcggt gatgcgagct cttccgaagg t 41

<210> SEQ ID NO 29

<211> LENGTH: 41

<212> TYPE: DNA

<213> ORGANISM: Artificial Sequence

<220> FEATURE:

<223> OTHER INFORMATION: Description of Artificial Sequence: Synthetic primer

<400> SEQUENCE: 29

cctctctatg ggcagtcggt gatgttgacc ctgtgggcag a 41

<210> SEQ ID NO 30

<211> LENGTH: 43

<212> TYPE: DNA

<213> ORGANISM: Artificial Sequence

<220> FEATURE:

<223> OTHER INFORMATION: Description of Artificial Sequence: Synthetic primer

<400> SEQUENCE: 30

cctctctatg ggcagtcggt gattcctctg ggtagtgca g ctt 43

<210> SEQ ID NO 31

<211> LENGTH: 43

<212> TYPE: DNA

<213> ORGANISM: Artificial Sequence

<220> FEATURE:

<223> OTHER INFORMATION: Description of Artificial Sequence: Synthetic primer

<400> SEQUENCE: 31

cctctctatg ggcagtcggt gatcagagcc gaattcagaa agc 43

<210> SEQ ID NO 32

<211> LENGTH: 41

<212> TYPE: DNA

<213> ORGANISM: Artificial Sequence

<220> FEATURE:

<223> OTHER INFORMATION: Description of Artificial Sequence: Synthetic primer

<400> SEQUENCE: 32

cctctctatg ggcagtcggt gatggggtgg tgtccaggag t 41

<210> SEQ ID NO 33

<211> LENGTH: 41

<212> TYPE: DNA

<213> ORGANISM: Artificial Sequence

<220> FEATURE:

-continued

<223> OTHER INFORMATION: Description of Artificial Sequence: Synthetic primer

<400> SEQUENCE: 33

cctctctatg ggcagtcggt gatcattttg aggcgcgaga a 41

<210> SEQ ID NO 34

<211> LENGTH: 41

<212> TYPE: DNA

<213> ORGANISM: Artificial Sequence

<220> FEATURE:

<223> OTHER INFORMATION: Description of Artificial Sequence: Synthetic primer

<400> SEQUENCE: 34

cctctctatg ggcagtcggt gatgaggcag ggagcagagg t 41

<210> SEQ ID NO 35

<211> LENGTH: 42

<212> TYPE: DNA

<213> ORGANISM: Artificial Sequence

<220> FEATURE:

<223> OTHER INFORMATION: Description of Artificial Sequence: Synthetic primer

<400> SEQUENCE: 35

cctctctatg ggcagtcggt gatcctttaa atcggggagg ag 42

<210> SEQ ID NO 36

<211> LENGTH: 44

<212> TYPE: DNA

<213> ORGANISM: Artificial Sequence

<220> FEATURE:

<223> OTHER INFORMATION: Description of Artificial Sequence: Synthetic primer

<400> SEQUENCE: 36

cctctctatg ggcagtcggt gatgagcagg aagggcaggt ataa 44

<210> SEQ ID NO 37

<211> LENGTH: 42

<212> TYPE: DNA

<213> ORGANISM: Artificial Sequence

<220> FEATURE:

<223> OTHER INFORMATION: Description of Artificial Sequence: Synthetic primer

<400> SEQUENCE: 37

cctctctatg ggcagtcggt gatctggcac aggtgactc at 42

<210> SEQ ID NO 38

<211> LENGTH: 41

<212> TYPE: DNA

<213> ORGANISM: Artificial Sequence

<220> FEATURE:

<223> OTHER INFORMATION: Description of Artificial Sequence: Synthetic primer

<400> SEQUENCE: 38

cctctctatg ggcagtcggt gatggatgcc tccagcagtg a 41

<210> SEQ ID NO 39

<211> LENGTH: 45

<212> TYPE: DNA

-continued

<213> ORGANISM: Artificial Sequence
<220> FEATURE:
<223> OTHER INFORMATION: Description of Artificial Sequence: Synthetic primer

<400> SEQUENCE: 39

cctctctatg ggcagtcggt gatagcctca taagggatgt actcg 45

<210> SEQ ID NO 40
<211> LENGTH: 43
<212> TYPE: DNA
<213> ORGANISM: Artificial Sequence
<220> FEATURE:
<223> OTHER INFORMATION: Description of Artificial Sequence: Synthetic primer

<400> SEQUENCE: 40

cctctctatg ggcagtcggt gatccttgaa aggactgcct gag 43

<210> SEQ ID NO 41
<211> LENGTH: 43
<212> TYPE: DNA
<213> ORGANISM: Artificial Sequence
<220> FEATURE:
<223> OTHER INFORMATION: Description of Artificial Sequence: Synthetic primer

<400> SEQUENCE: 41

cctctctatg ggcagtcggt gatctcctca cctggcactt gat 43

<210> SEQ ID NO 42
<211> LENGTH: 43
<212> TYPE: DNA
<213> ORGANISM: Artificial Sequence
<220> FEATURE:
<223> OTHER INFORMATION: Description of Artificial Sequence: Synthetic primer

<400> SEQUENCE: 42

cctctctatg ggcagtcggt gatagagccc acacaggaga gag 43

<210> SEQ ID NO 43
<211> LENGTH: 42
<212> TYPE: DNA
<213> ORGANISM: Artificial Sequence
<220> FEATURE:
<223> OTHER INFORMATION: Description of Artificial Sequence: Synthetic primer

<400> SEQUENCE: 43

cctctctatg ggcagtcggt gatatcagcg agtccacact gg 42

<210> SEQ ID NO 44
<211> LENGTH: 47
<212> TYPE: DNA
<213> ORGANISM: Artificial Sequence
<220> FEATURE:
<223> OTHER INFORMATION: Description of Artificial Sequence: Synthetic primer

<400> SEQUENCE: 44

cctctctatg ggcagtcggt gatttgaatc ataagagact ccattgc 47

<210> SEQ ID NO 45

-continued

<211> LENGTH: 58
<212> TYPE: DNA
<213> ORGANISM: Artificial Sequence
<220> FEATURE:
<223> OTHER INFORMATION: Description of Artificial Sequence: Synthetic primer

<400> SEQUENCE: 45

ccatctcatc cctgcgtgtc tccgactcag ttaacggacg cgccagatgg gtaagtgc 58

<210> SEQ ID NO 46
<211> LENGTH: 59
<212> TYPE: DNA
<213> ORGANISM: Artificial Sequence
<220> FEATURE:
<223> OTHER INFORMATION: Description of Artificial Sequence: Synthetic primer

<400> SEQUENCE: 46

ccatctcatc cctgcgtgtc tccgactcag tccggcttac tgcaaatcag tggctctcc 59

<210> SEQ ID NO 47
<211> LENGTH: 58
<212> TYPE: DNA
<213> ORGANISM: Artificial Sequence
<220> FEATURE:
<223> OTHER INFORMATION: Description of Artificial Sequence: Synthetic primer

<400> SEQUENCE: 47

ccatctcatc cctgcgtgtc tccgactcag tctcattcag ctcccttctt gccacctg 58

<210> SEQ ID NO 48
<211> LENGTH: 65
<212> TYPE: DNA
<213> ORGANISM: Artificial Sequence
<220> FEATURE:
<223> OTHER INFORMATION: Description of Artificial Sequence: Synthetic primer

<400> SEQUENCE: 48

ccatctcatc cctgcgtgtc tccgactcag gcggtcatac acatgtgata cttttgggaa 60
tgaag 65

<210> SEQ ID NO 49
<211> LENGTH: 62
<212> TYPE: DNA
<213> ORGANISM: Artificial Sequence
<220> FEATURE:
<223> OTHER INFORMATION: Description of Artificial Sequence: Synthetic primer

<400> SEQUENCE: 49

ccatctcatc cctgcgtgtc tccgactcag taggacgttc cttctgaaca ccaaattgga 60
aa 62

<210> SEQ ID NO 50
<211> LENGTH: 61
<212> TYPE: DNA
<213> ORGANISM: Artificial Sequence
<220> FEATURE:
<223> OTHER INFORMATION: Description of Artificial Sequence: Synthetic primer

-continued

<400> SEQUENCE: 50

ccatctcattc cctgcgtgtc tccgactcag acgacgcaac tgtaagagc ccagagggtc 60

a 61

<210> SEQ ID NO 51

<211> LENGTH: 58

<212> TYPE: DNA

<213> ORGANISM: Artificial Sequence

<220> FEATURE:

<223> OTHER INFORMATION: Description of Artificial Sequence: Synthetic primer

<400> SEQUENCE: 51

ccatctcattc cctgcgtgtc tccgactcag cttctcggac ggggcccta ctcttga 58

<210> SEQ ID NO 52

<211> LENGTH: 58

<212> TYPE: DNA

<213> ORGANISM: Artificial Sequence

<220> FEATURE:

<223> OTHER INFORMATION: Description of Artificial Sequence: Synthetic primer

<400> SEQUENCE: 52

ccatctcattc cctgcgtgtc tccgactcag cattgccgtt ctgcagccc ctacactg 58

<210> SEQ ID NO 53

<211> LENGTH: 58

<212> TYPE: DNA

<213> ORGANISM: Artificial Sequence

<220> FEATURE:

<223> OTHER INFORMATION: Description of Artificial Sequence: Synthetic primer

<400> SEQUENCE: 53

ccatctcattc cctgcgtgtc tccgactcag cgagccagaa tgactgccct ctgctgtg 58

<210> SEQ ID NO 54

<211> LENGTH: 60

<212> TYPE: DNA

<213> ORGANISM: Artificial Sequence

<220> FEATURE:

<223> OTHER INFORMATION: Description of Artificial Sequence: Synthetic primer

<400> SEQUENCE: 54

ccatctcattc cctgcgtgtc tccgactcag tgaggacggc agtgacgatg agcaggaggt 60

<210> SEQ ID NO 55

<211> LENGTH: 58

<212> TYPE: DNA

<213> ORGANISM: Artificial Sequence

<220> FEATURE:

<223> OTHER INFORMATION: Description of Artificial Sequence: Synthetic primer

<400> SEQUENCE: 55

ccatctcattc cctgcgtgtc tccgactcag gcctgcgcag gccaatccc attgacca 58

<210> SEQ ID NO 56

<211> LENGTH: 60

<212> TYPE: DNA

<213> ORGANISM: Artificial Sequence

-continued

<220> FEATURE:
<223> OTHER INFORMATION: Description of Artificial Sequence: Synthetic primer

<400> SEQUENCE: 56

ccatctcatc cctgcgtgtc tccgactcag gttgacgtct ccaagcttcc tgaaccagac 60

<210> SEQ ID NO 57
<211> LENGTH: 58
<212> TYPE: DNA
<213> ORGANISM: Artificial Sequence
<220> FEATURE:
<223> OTHER INFORMATION: Description of Artificial Sequence: Synthetic primer

<400> SEQUENCE: 57

ccatctcatc cctgcgtgtc tccgactcag gagatcgatt ctagtggggg cattccaa 58

<210> SEQ ID NO 58
<211> LENGTH: 58
<212> TYPE: DNA
<213> ORGANISM: Artificial Sequence
<220> FEATURE:
<223> OTHER INFORMATION: Description of Artificial Sequence: Synthetic primer

<400> SEQUENCE: 58

ccatctcatc cctgcgtgtc tccgactcag agttcgagcc ctctaggggc cgtttcct 58

<210> SEQ ID NO 59
<211> LENGTH: 59
<212> TYPE: DNA
<213> ORGANISM: Artificial Sequence
<220> FEATURE:
<223> OTHER INFORMATION: Description of Artificial Sequence: Synthetic primer

<400> SEQUENCE: 59

ccatctcatc cctgcgtgtc tccgactcag ctcaggetca tcagccttcc ctcgtctta 59

<210> SEQ ID NO 60
<211> LENGTH: 60
<212> TYPE: DNA
<213> ORGANISM: Artificial Sequence
<220> FEATURE:
<223> OTHER INFORMATION: Description of Artificial Sequence: Synthetic primer

<400> SEQUENCE: 60

ccatctcatc cctgcgtgtc tccgactcag ggcaatataa tccacataac tcgcttgcag 60

<210> SEQ ID NO 61
<211> LENGTH: 61
<212> TYPE: DNA
<213> ORGANISM: Artificial Sequence
<220> FEATURE:
<223> OTHER INFORMATION: Description of Artificial Sequence: Synthetic primer

<400> SEQUENCE: 61

ccatctcatc cctgcgtgtc tccgactcag ggtactcatg gaactgtagc ccagacactg 60

c 61

-continued

<210> SEQ ID NO 62
<211> LENGTH: 63
<212> TYPE: DNA
<213> ORGANISM: Artificial Sequence
<220> FEATURE:
<223> OTHER INFORMATION: Description of Artificial Sequence: Synthetic primer

<400> SEQUENCE: 62

ccatctcattc cctgcgtgtc tccgactcag tctggttcaa acaaagctgg aaactcttcc 60

cta 63

<210> SEQ ID NO 63
<211> LENGTH: 60
<212> TYPE: DNA
<213> ORGANISM: Artificial Sequence
<220> FEATURE:
<223> OTHER INFORMATION: Description of Artificial Sequence: Synthetic primer

<400> SEQUENCE: 63

ccatctcattc cctgcgtgtc tccgactcag gtcctataag ccaacaagcc caacaagtcc 60

<210> SEQ ID NO 64
<211> LENGTH: 60
<212> TYPE: DNA
<213> ORGANISM: Artificial Sequence
<220> FEATURE:
<223> OTHER INFORMATION: Description of Artificial Sequence: Synthetic primer

<400> SEQUENCE: 64

ccatctcattc cctgcgtgtc tccgactcag gtcagcctcc gaatgaccgg ctgtctgttt 60

<210> SEQ ID NO 65
<211> LENGTH: 59
<212> TYPE: DNA
<213> ORGANISM: Artificial Sequence
<220> FEATURE:
<223> OTHER INFORMATION: Description of Artificial Sequence: Synthetic primer

<400> SEQUENCE: 65

ccatctcattc cctgcgtgtc tccgactcag ttcaagctcg aaagtggcac caccaacaa 59

<210> SEQ ID NO 66
<211> LENGTH: 60
<212> TYPE: DNA
<213> ORGANISM: Artificial Sequence
<220> FEATURE:
<223> OTHER INFORMATION: Description of Artificial Sequence: Synthetic primer

<400> SEQUENCE: 66

ccatctcattc cctgcgtgtc tccgactcag gtaccagcgc cttgtagcgc ttcccacagt 60

<210> SEQ ID NO 67
<211> LENGTH: 61
<212> TYPE: DNA
<213> ORGANISM: Artificial Sequence
<220> FEATURE:
<223> OTHER INFORMATION: Description of Artificial Sequence: Synthetic primer

<400> SEQUENCE: 67

-continued

ccatctcadc cctgcgtgtc tccgactcag tectattcgg agcttcttcc cacaatcctc 60

a 61

<210> SEQ ID NO 68

<211> LENGTH: 67

<212> TYPE: DNA

<213> ORGANISM: Artificial Sequence

<220> FEATURE:

<223> OTHER INFORMATION: Description of Artificial Sequence: Synthetic primer

<400> SEQUENCE: 68

ccatctcadc cctgcgtgtc tccgactcag gccagcgatt ctgtacccca taaatatgta 60

caacact 67

<210> SEQ ID NO 69

<211> LENGTH: 58

<212> TYPE: DNA

<213> ORGANISM: Artificial Sequence

<220> FEATURE:

<223> OTHER INFORMATION: Description of Artificial Sequence: Synthetic primer

<400> SEQUENCE: 69

ccatctcadc cctgcgtgtc tccgactcag acctagactg cgccagatgg gtaagtgc 58

<210> SEQ ID NO 70

<211> LENGTH: 59

<212> TYPE: DNA

<213> ORGANISM: Artificial Sequence

<220> FEATURE:

<223> OTHER INFORMATION: Description of Artificial Sequence: Synthetic primer

<400> SEQUENCE: 70

ccatctcadc cctgcgtgtc tccgactcag actgggttcgc tgcaaatcag tggctctcc 59

<210> SEQ ID NO 71

<211> LENGTH: 58

<212> TYPE: DNA

<213> ORGANISM: Artificial Sequence

<220> FEATURE:

<223> OTHER INFORMATION: Description of Artificial Sequence: Synthetic primer

<400> SEQUENCE: 71

ccatctcadc cctgcgtgtc tccgactcag ccatattagg ctcccttctc gccacctg 58

<210> SEQ ID NO 72

<211> LENGTH: 65

<212> TYPE: DNA

<213> ORGANISM: Artificial Sequence

<220> FEATURE:

<223> OTHER INFORMATION: Description of Artificial Sequence: Synthetic primer

<400> SEQUENCE: 72

ccatctcadc cctgcgtgtc tccgactcag gctcgtcagc acatgtgata cttttgggaa 60

tgaag 65

<210> SEQ ID NO 73

-continued

<211> LENGTH: 62
<212> TYPE: DNA
<213> ORGANISM: Artificial Sequence
<220> FEATURE:
<223> OTHER INFORMATION: Description of Artificial Sequence: Synthetic primer

<400> SEQUENCE: 73

ccatctcatc cctgcgtgtc tccgactcag cgtaatgacg cttctgaaca ccaaattgga 60
aa 62

<210> SEQ ID NO 74
<211> LENGTH: 61
<212> TYPE: DNA
<213> ORGANISM: Artificial Sequence
<220> FEATURE:
<223> OTHER INFORMATION: Description of Artificial Sequence: Synthetic primer

<400> SEQUENCE: 74

ccatctcatc cctgcgtgtc tccgactcag cggcgctga tgttaagagc ccagagggtc 60
a 61

<210> SEQ ID NO 75
<211> LENGTH: 58
<212> TYPE: DNA
<213> ORGANISM: Artificial Sequence
<220> FEATURE:
<223> OTHER INFORMATION: Description of Artificial Sequence: Synthetic primer

<400> SEQUENCE: 75

ccatctcatc cctgcgtgtc tccgactcag cgcaagata ggggcccta ctcttga 58

<210> SEQ ID NO 76
<211> LENGTH: 58
<212> TYPE: DNA
<213> ORGANISM: Artificial Sequence
<220> FEATURE:
<223> OTHER INFORMATION: Description of Artificial Sequence: Synthetic primer

<400> SEQUENCE: 76

ccatctcatc cctgcgtgtc tccgactcag gaaccgcaga ctgcagccc ctacactg 58

<210> SEQ ID NO 77
<211> LENGTH: 58
<212> TYPE: DNA
<213> ORGANISM: Artificial Sequence
<220> FEATURE:
<223> OTHER INFORMATION: Description of Artificial Sequence: Synthetic primer

<400> SEQUENCE: 77

ccatctcatc cctgcgtgtc tccgactcag ttggcagaga tgactgccct ctgctgtg 58

<210> SEQ ID NO 78
<211> LENGTH: 60
<212> TYPE: DNA
<213> ORGANISM: Artificial Sequence
<220> FEATURE:
<223> OTHER INFORMATION: Description of Artificial Sequence: Synthetic primer

-continued

<400> SEQUENCE: 78

ccatctcattc cctgcgtgtc tccgactcag gcattctctgc agtgacgatg agcaggaggt 60

<210> SEQ ID NO 79

<211> LENGTH: 58

<212> TYPE: DNA

<213> ORGANISM: Artificial Sequence

<220> FEATURE:

<223> OTHER INFORMATION: Description of Artificial Sequence: Synthetic primer

<400> SEQUENCE: 79

ccatctcattc cctgcgtgtc tccgactcag ttggaccgca gccaatcccc attgacca 58

<210> SEQ ID NO 80

<211> LENGTH: 60

<212> TYPE: DNA

<213> ORGANISM: Artificial Sequence

<220> FEATURE:

<223> OTHER INFORMATION: Description of Artificial Sequence: Synthetic primer

<400> SEQUENCE: 80

ccatctcattc cctgcgtgtc tccgactcag gcagaacgtc ccaagcttcc tgaaccagac 60

<210> SEQ ID NO 81

<211> LENGTH: 58

<212> TYPE: DNA

<213> ORGANISM: Artificial Sequence

<220> FEATURE:

<223> OTHER INFORMATION: Description of Artificial Sequence: Synthetic primer

<400> SEQUENCE: 81

ccatctcattc cctgcgtgtc tccgactcag aacttcgagc ctagtggggg cattccaa 58

<210> SEQ ID NO 82

<211> LENGTH: 58

<212> TYPE: DNA

<213> ORGANISM: Artificial Sequence

<220> FEATURE:

<223> OTHER INFORMATION: Description of Artificial Sequence: Synthetic primer

<400> SEQUENCE: 82

ccatctcattc cctgcgtgtc tccgactcag gctcctagag ctctaggggc cgtttcct 58

<210> SEQ ID NO 83

<211> LENGTH: 59

<212> TYPE: DNA

<213> ORGANISM: Artificial Sequence

<220> FEATURE:

<223> OTHER INFORMATION: Description of Artificial Sequence: Synthetic primer

<400> SEQUENCE: 83

ccatctcattc cctgcgtgtc tccgactcag tatctagctt tcagccttcc ctgcctcta 59

<210> SEQ ID NO 84

<211> LENGTH: 60

<212> TYPE: DNA

<213> ORGANISM: Artificial Sequence

<220> FEATURE:

<223> OTHER INFORMATION: Description of Artificial Sequence: Synthetic

-continued

primer

<400> SEQUENCE: 84

ccatctcatc cctgcgtgtc tccgactcag gagtattggc tccacataac tcgcttgacg 60

<210> SEQ ID NO 85
<211> LENGTH: 61
<212> TYPE: DNA
<213> ORGANISM: Artificial Sequence
<220> FEATURE:
<223> OTHER INFORMATION: Description of Artificial Sequence: Synthetic primer

<400> SEQUENCE: 85

ccatctcatc cctgcgtgtc tccgactcag cctgagctca gaactgtacg ccagacactg 60

c 61

<210> SEQ ID NO 86
<211> LENGTH: 63
<212> TYPE: DNA
<213> ORGANISM: Artificial Sequence
<220> FEATURE:
<223> OTHER INFORMATION: Description of Artificial Sequence: Synthetic primer

<400> SEQUENCE: 86

ccatctcatc cctgcgtgtc tccgactcag caggcgagta acaaagctgg aaactcttcc 60

cta 63

<210> SEQ ID NO 87
<211> LENGTH: 60
<212> TYPE: DNA
<213> ORGANISM: Artificial Sequence
<220> FEATURE:
<223> OTHER INFORMATION: Description of Artificial Sequence: Synthetic primer

<400> SEQUENCE: 87

ccatctcatc cctgcgtgtc tccgactcag gcaggcagag ccaacaagcc caacaagttc 60

<210> SEQ ID NO 88
<211> LENGTH: 60
<212> TYPE: DNA
<213> ORGANISM: Artificial Sequence
<220> FEATURE:
<223> OTHER INFORMATION: Description of Artificial Sequence: Synthetic primer

<400> SEQUENCE: 88

ccatctcatc cctgcgtgtc tccgactcag gctcgatac gaatgaccgg ctgtctgttt 60

<210> SEQ ID NO 89
<211> LENGTH: 59
<212> TYPE: DNA
<213> ORGANISM: Artificial Sequence
<220> FEATURE:
<223> OTHER INFORMATION: Description of Artificial Sequence: Synthetic primer

<400> SEQUENCE: 89

ccatctcatc cctgcgtgtc tccgactcag cgatgattat aaagtggcac caccaacaa 59

-continued

<210> SEQ ID NO 90
<211> LENGTH: 60
<212> TYPE: DNA
<213> ORGANISM: Artificial Sequence
<220> FEATURE:
<223> OTHER INFORMATION: Description of Artificial Sequence: Synthetic primer

<400> SEQUENCE: 90

ccatctcatt cctgcgtgtc tccgactcag gacggctggc ctgttagcgc ttcccacagt 60

<210> SEQ ID NO 91
<211> LENGTH: 61
<212> TYPE: DNA
<213> ORGANISM: Artificial Sequence
<220> FEATURE:
<223> OTHER INFORMATION: Description of Artificial Sequence: Synthetic primer

<400> SEQUENCE: 91

ccatctcatt cctgcgtgtc tccgactcag ggagcctgag agcttctttc cacaatcctc 60

a 61

<210> SEQ ID NO 92
<211> LENGTH: 67
<212> TYPE: DNA
<213> ORGANISM: Artificial Sequence
<220> FEATURE:
<223> OTHER INFORMATION: Description of Artificial Sequence: Synthetic primer

<400> SEQUENCE: 92

ccatctcatt cctgcgtgtc tccgactcag cctgactgct ctgtacccca taaatatgta 60

caacact 67

<210> SEQ ID NO 93
<211> LENGTH: 58
<212> TYPE: DNA
<213> ORGANISM: Artificial Sequence
<220> FEATURE:
<223> OTHER INFORMATION: Description of Artificial Sequence: Synthetic primer

<400> SEQUENCE: 93

ccatctcatt cctgcgtgtc tccgactcag acggctgacg cgccagatgg gtaagtgc 58

<210> SEQ ID NO 94
<211> LENGTH: 59
<212> TYPE: DNA
<213> ORGANISM: Artificial Sequence
<220> FEATURE:
<223> OTHER INFORMATION: Description of Artificial Sequence: Synthetic primer

<400> SEQUENCE: 94

ccatctcatt cctgcgtgtc tccgactcag taacctatgc tgcaaatcag tggctctcc 59

<210> SEQ ID NO 95
<211> LENGTH: 58
<212> TYPE: DNA
<213> ORGANISM: Artificial Sequence
<220> FEATURE:
<223> OTHER INFORMATION: Description of Artificial Sequence: Synthetic primer

-continued

<400> SEQUENCE: 95

ccatctcatc cctgcgtgtc tccgactcag tcttgccctc ctcccttctt gccacctg 58

<210> SEQ ID NO 96

<211> LENGTH: 65

<212> TYPE: DNA

<213> ORGANISM: Artificial Sequence

<220> FEATURE:

<223> OTHER INFORMATION: Description of Artificial Sequence: Synthetic primer

<400> SEQUENCE: 96

ccatctcatc cctgcgtgtc tccgactcag ttcttagatt acatgtgata cttttgggaa 60

tgaag 65

<210> SEQ ID NO 97

<211> LENGTH: 62

<212> TYPE: DNA

<213> ORGANISM: Artificial Sequence

<220> FEATURE:

<223> OTHER INFORMATION: Description of Artificial Sequence: Synthetic primer

<400> SEQUENCE: 97

ccatctcatc cctgcgtgtc tccgactcag tcattctcatt cttctgaaca ccaaattgga 60

aa 62

<210> SEQ ID NO 98

<211> LENGTH: 61

<212> TYPE: DNA

<213> ORGANISM: Artificial Sequence

<220> FEATURE:

<223> OTHER INFORMATION: Description of Artificial Sequence: Synthetic primer

<400> SEQUENCE: 98

ccatctcatc cctgcgtgtc tccgactcag tctccgctcg tgttaagagc ccagagggtc 60

a 61

<210> SEQ ID NO 99

<211> LENGTH: 58

<212> TYPE: DNA

<213> ORGANISM: Artificial Sequence

<220> FEATURE:

<223> OTHER INFORMATION: Description of Artificial Sequence: Synthetic primer

<400> SEQUENCE: 99

ccatctcatc cctgcgtgtc tccgactcag tgccatatgc ggggccccta ctctttga 58

<210> SEQ ID NO 100

<211> LENGTH: 58

<212> TYPE: DNA

<213> ORGANISM: Artificial Sequence

<220> FEATURE:

<223> OTHER INFORMATION: Description of Artificial Sequence: Synthetic primer

<400> SEQUENCE: 100

ccatctcatc cctgcgtgtc tccgactcag taaggcctct ctgcagccc ctacactg 58

-continued

<210> SEQ ID NO 101
<211> LENGTH: 58
<212> TYPE: DNA
<213> ORGANISM: Artificial Sequence
<220> FEATURE:
<223> OTHER INFORMATION: Description of Artificial Sequence: Synthetic primer

<400> SEQUENCE: 101

ccatctcatc cctgcgtgtc tccgactcag gagtaggccg tgactgccct ctgctgtg 58

<210> SEQ ID NO 102
<211> LENGTH: 60
<212> TYPE: DNA
<213> ORGANISM: Artificial Sequence
<220> FEATURE:
<223> OTHER INFORMATION: Description of Artificial Sequence: Synthetic primer

<400> SEQUENCE: 102

ccatctcatc cctgcgtgtc tccgactcag gcaataagct agtgacgatg agcaggaggt 60

<210> SEQ ID NO 103
<211> LENGTH: 58
<212> TYPE: DNA
<213> ORGANISM: Artificial Sequence
<220> FEATURE:
<223> OTHER INFORMATION: Description of Artificial Sequence: Synthetic primer

<400> SEQUENCE: 103

ccatctcatc cctgcgtgtc tccgactcag ggcgttgcaa gccaatccc attgacca 58

<210> SEQ ID NO 104
<211> LENGTH: 60
<212> TYPE: DNA
<213> ORGANISM: Artificial Sequence
<220> FEATURE:
<223> OTHER INFORMATION: Description of Artificial Sequence: Synthetic primer

<400> SEQUENCE: 104

ccatctcatc cctgcgtgtc tccgactcag ccaagaagcg ccaagcttcc tgaaccagac 60

<210> SEQ ID NO 105
<211> LENGTH: 58
<212> TYPE: DNA
<213> ORGANISM: Artificial Sequence
<220> FEATURE:
<223> OTHER INFORMATION: Description of Artificial Sequence: Synthetic primer

<400> SEQUENCE: 105

ccatctcatc cctgcgtgtc tccgactcag gggtacctcg ctagtggggg cattccaa 58

<210> SEQ ID NO 106
<211> LENGTH: 58
<212> TYPE: DNA
<213> ORGANISM: Artificial Sequence
<220> FEATURE:
<223> OTHER INFORMATION: Description of Artificial Sequence: Synthetic primer

<400> SEQUENCE: 106

-continued

ccatctcatc cctgcgtgtc tccgactcag ctccgcctta ctctagggcg cgtttcct 58

<210> SEQ ID NO 107
<211> LENGTH: 59
<212> TYPE: DNA
<213> ORGANISM: Artificial Sequence
<220> FEATURE:
<223> OTHER INFORMATION: Description of Artificial Sequence: Synthetic primer

<400> SEQUENCE: 107

ccatctcatc cctgcgtgtc tccgactcag ctccagagat tcagcctttc ctgcgtcta 59

<210> SEQ ID NO 108
<211> LENGTH: 60
<212> TYPE: DNA
<213> ORGANISM: Artificial Sequence
<220> FEATURE:
<223> OTHER INFORMATION: Description of Artificial Sequence: Synthetic primer

<400> SEQUENCE: 108

ccatctcatc cctgcgtgtc tccgactcag gtcgaggtag tccacataac tcgcttgca 60

<210> SEQ ID NO 109
<211> LENGTH: 61
<212> TYPE: DNA
<213> ORGANISM: Artificial Sequence
<220> FEATURE:
<223> OTHER INFORMATION: Description of Artificial Sequence: Synthetic primer

<400> SEQUENCE: 109

ccatctcatc cctgcgtgtc tccgactcag tatggacctg gaactgtagc ccagacactg 60

c 61

<210> SEQ ID NO 110
<211> LENGTH: 63
<212> TYPE: DNA
<213> ORGANISM: Artificial Sequence
<220> FEATURE:
<223> OTHER INFORMATION: Description of Artificial Sequence: Synthetic primer

<400> SEQUENCE: 110

ccatctcatc cctgcgtgtc tccgactcag tacctgctag acaaagctgg aaactcttcc 60

cta 63

<210> SEQ ID NO 111
<211> LENGTH: 60
<212> TYPE: DNA
<213> ORGANISM: Artificial Sequence
<220> FEATURE:
<223> OTHER INFORMATION: Description of Artificial Sequence: Synthetic primer

<400> SEQUENCE: 111

ccatctcatc cctgcgtgtc tccgactcag ccgcgaccga ccaacaagcc caacaagttc 60

<210> SEQ ID NO 112
<211> LENGTH: 60
<212> TYPE: DNA
<213> ORGANISM: Artificial Sequence

-continued

<220> FEATURE:
<223> OTHER INFORMATION: Description of Artificial Sequence: Synthetic primer

<400> SEQUENCE: 112

ccatctcattc cctgcgtgtc tccgactcag gttgaacgtt gaatgaccgg ctgtctgttt 60

<210> SEQ ID NO 113
<211> LENGTH: 59
<212> TYPE: DNA
<213> ORGANISM: Artificial Sequence
<220> FEATURE:
<223> OTHER INFORMATION: Description of Artificial Sequence: Synthetic primer

<400> SEQUENCE: 113

ccatctcattc cctgcgtgtc tccgactcag tgccaacgca aaagtggcac caccaacaa 59

<210> SEQ ID NO 114
<211> LENGTH: 60
<212> TYPE: DNA
<213> ORGANISM: Artificial Sequence
<220> FEATURE:
<223> OTHER INFORMATION: Description of Artificial Sequence: Synthetic primer

<400> SEQUENCE: 114

ccatctcattc cctgcgtgtc tccgactcag ggattgacct ctgttagcgc ttcccacagt 60

<210> SEQ ID NO 115
<211> LENGTH: 61
<212> TYPE: DNA
<213> ORGANISM: Artificial Sequence
<220> FEATURE:
<223> OTHER INFORMATION: Description of Artificial Sequence: Synthetic primer

<400> SEQUENCE: 115

ccatctcattc cctgcgtgtc tccgactcag ggacggattc agcttctttc cacaatcctc 60

a 61

<210> SEQ ID NO 116
<211> LENGTH: 67
<212> TYPE: DNA
<213> ORGANISM: Artificial Sequence
<220> FEATURE:
<223> OTHER INFORMATION: Description of Artificial Sequence: Synthetic primer

<400> SEQUENCE: 116

ccatctcattc cctgcgtgtc tccgactcag tcttcgctcg ctgtaccocca taaatatgta 60

caacact 67

<210> SEQ ID NO 117
<211> LENGTH: 58
<212> TYPE: DNA
<213> ORGANISM: Artificial Sequence
<220> FEATURE:
<223> OTHER INFORMATION: Description of Artificial Sequence: Synthetic primer

<400> SEQUENCE: 117

ccatctcattc cctgcgtgtc tccgactcag agttcatggt cgccagatgg gtaagtgc 58

-continued

<210> SEQ ID NO 118
<211> LENGTH: 59
<212> TYPE: DNA
<213> ORGANISM: Artificial Sequence
<220> FEATURE:
<223> OTHER INFORMATION: Description of Artificial Sequence: Synthetic
primer

<400> SEQUENCE: 118

ccatctcatc cctgcgtgtc tccgactcag tatccattcc tgcaaatcag tggtctctcc 59

<210> SEQ ID NO 119
<211> LENGTH: 58
<212> TYPE: DNA
<213> ORGANISM: Artificial Sequence
<220> FEATURE:
<223> OTHER INFORMATION: Description of Artificial Sequence: Synthetic
primer

<400> SEQUENCE: 119

ccatctcatc cctgcgtgtc tccgactcag ggagagcgcg ctcccttctc gccacctg 58

<210> SEQ ID NO 120
<211> LENGTH: 65
<212> TYPE: DNA
<213> ORGANISM: Artificial Sequence
<220> FEATURE:
<223> OTHER INFORMATION: Description of Artificial Sequence: Synthetic
primer

<400> SEQUENCE: 120

ccatctcatc cctgcgtgtc tccgactcag cggaccttgg acatgtgata cttttgggaa 60
tgaag 65

<210> SEQ ID NO 121
<211> LENGTH: 62
<212> TYPE: DNA
<213> ORGANISM: Artificial Sequence
<220> FEATURE:
<223> OTHER INFORMATION: Description of Artificial Sequence: Synthetic
primer

<400> SEQUENCE: 121

ccatctcatc cctgcgtgtc tccgactcag ggcaatctcc cttctgaaca ccaaattgga 60
aa 62

<210> SEQ ID NO 122
<211> LENGTH: 61
<212> TYPE: DNA
<213> ORGANISM: Artificial Sequence
<220> FEATURE:
<223> OTHER INFORMATION: Description of Artificial Sequence: Synthetic
primer

<400> SEQUENCE: 122

ccatctcatc cctgcgtgtc tccgactcag aggattgatt tgtaagagc ccagaggttc 60
a 61

<210> SEQ ID NO 123
<211> LENGTH: 58
<212> TYPE: DNA

-continued

<213> ORGANISM: Artificial Sequence
<220> FEATURE:
<223> OTHER INFORMATION: Description of Artificial Sequence: Synthetic primer

<400> SEQUENCE: 123

ccatctcatc cctgcgtgtc tccgactcag gccgttgccct ggggcccta ctctttga 58

<210> SEQ ID NO 124
<211> LENGTH: 58
<212> TYPE: DNA
<213> ORGANISM: Artificial Sequence
<220> FEATURE:
<223> OTHER INFORMATION: Description of Artificial Sequence: Synthetic primer

<400> SEQUENCE: 124

ccatctcatc cctgcgtgtc tccgactcag aagtaagtcg ctgcagccc ctacactg 58

<210> SEQ ID NO 125
<211> LENGTH: 58
<212> TYPE: DNA
<213> ORGANISM: Artificial Sequence
<220> FEATURE:
<223> OTHER INFORMATION: Description of Artificial Sequence: Synthetic primer

<400> SEQUENCE: 125

ccatctcatc cctgcgtgtc tccgactcag tggcttaagg tgactgccct ctgctgtg 58

<210> SEQ ID NO 126
<211> LENGTH: 60
<212> TYPE: DNA
<213> ORGANISM: Artificial Sequence
<220> FEATURE:
<223> OTHER INFORMATION: Description of Artificial Sequence: Synthetic primer

<400> SEQUENCE: 126

ccatctcatc cctgcgtgtc tccgactcag ctcttcaga agtgacgatg agcaggaggt 60

<210> SEQ ID NO 127
<211> LENGTH: 58
<212> TYPE: DNA
<213> ORGANISM: Artificial Sequence
<220> FEATURE:
<223> OTHER INFORMATION: Description of Artificial Sequence: Synthetic primer

<400> SEQUENCE: 127

ccatctcatc cctgcgtgtc tccgactcag cgttcttcaa gccaatccc attgacca 58

<210> SEQ ID NO 128
<211> LENGTH: 60
<212> TYPE: DNA
<213> ORGANISM: Artificial Sequence
<220> FEATURE:
<223> OTHER INFORMATION: Description of Artificial Sequence: Synthetic primer

<400> SEQUENCE: 128

ccatctcatc cctgcgtgtc tccgactcag caacggctgc ccaagcttcc tgaaccagac 60

<210> SEQ ID NO 129

-continued

<211> LENGTH: 58
<212> TYPE: DNA
<213> ORGANISM: Artificial Sequence
<220> FEATURE:
<223> OTHER INFORMATION: Description of Artificial Sequence: Synthetic primer
<400> SEQUENCE: 129
ccatctcatc cctgcgtgtc tccgactcag gcaagtaacc ctagtggggg cattccaa 58

<210> SEQ ID NO 130
<211> LENGTH: 58
<212> TYPE: DNA
<213> ORGANISM: Artificial Sequence
<220> FEATURE:
<223> OTHER INFORMATION: Description of Artificial Sequence: Synthetic primer
<400> SEQUENCE: 130
ccatctcatc cctgcgtgtc tccgactcag gtccatagtc ctctagggcg cgtttcct 58

<210> SEQ ID NO 131
<211> LENGTH: 59
<212> TYPE: DNA
<213> ORGANISM: Artificial Sequence
<220> FEATURE:
<223> OTHER INFORMATION: Description of Artificial Sequence: Synthetic primer
<400> SEQUENCE: 131
ccatctcatc cctgcgtgtc tccgactcag acggcgagcc tcagcctttc ctcgtctta 59

<210> SEQ ID NO 132
<211> LENGTH: 60
<212> TYPE: DNA
<213> ORGANISM: Artificial Sequence
<220> FEATURE:
<223> OTHER INFORMATION: Description of Artificial Sequence: Synthetic primer
<400> SEQUENCE: 132
ccatctcatc cctgcgtgtc tccgactcag gtatggtcgg tccacataac tcgcttcag 60

<210> SEQ ID NO 133
<211> LENGTH: 61
<212> TYPE: DNA
<213> ORGANISM: Artificial Sequence
<220> FEATURE:
<223> OTHER INFORMATION: Description of Artificial Sequence: Synthetic primer
<400> SEQUENCE: 133
ccatctcatc cctgcgtgtc tccgactcag tcggttatcc gaactgtagc ccagacactg 60
c 61

<210> SEQ ID NO 134
<211> LENGTH: 63
<212> TYPE: DNA
<213> ORGANISM: Artificial Sequence
<220> FEATURE:
<223> OTHER INFORMATION: Description of Artificial Sequence: Synthetic primer
<400> SEQUENCE: 134

-continued

ccatctcatc cctgcgtgtc tccgactcag gcggtcgata acaaagctgg aaactcttcc 60
cta 63

<210> SEQ ID NO 135
<211> LENGTH: 60
<212> TYPE: DNA
<213> ORGANISM: Artificial Sequence
<220> FEATURE:
<223> OTHER INFORMATION: Description of Artificial Sequence: Synthetic
primer

<400> SEQUENCE: 135

ccatctcatc cctgcgtgtc tccgactcag tcctcagtat ccaacaagcc caacaagttc 60

<210> SEQ ID NO 136
<211> LENGTH: 60
<212> TYPE: DNA
<213> ORGANISM: Artificial Sequence
<220> FEATURE:
<223> OTHER INFORMATION: Description of Artificial Sequence: Synthetic
primer

<400> SEQUENCE: 136

ccatctcatc cctgcgtgtc tccgactcag accgttctctg gaatgaccgg ctgtctgttt 60

<210> SEQ ID NO 137
<211> LENGTH: 59
<212> TYPE: DNA
<213> ORGANISM: Artificial Sequence
<220> FEATURE:
<223> OTHER INFORMATION: Description of Artificial Sequence: Synthetic
primer

<400> SEQUENCE: 137

ccatctcatc cctgcgtgtc tccgactcag gectgtcttt aaagtggcac caccaacaa 59

<210> SEQ ID NO 138
<211> LENGTH: 60
<212> TYPE: DNA
<213> ORGANISM: Artificial Sequence
<220> FEATURE:
<223> OTHER INFORMATION: Description of Artificial Sequence: Synthetic
primer

<400> SEQUENCE: 138

ccatctcatc cctgcgtgtc tccgactcag agcgtaacca cttgtagcgc ttcccacagt 60

<210> SEQ ID NO 139
<211> LENGTH: 61
<212> TYPE: DNA
<213> ORGANISM: Artificial Sequence
<220> FEATURE:
<223> OTHER INFORMATION: Description of Artificial Sequence: Synthetic
primer

<400> SEQUENCE: 139

ccatctcatc cctgcgtgtc tccgactcag ttgectgatg agcttctttc cacaatcctc 60

a 61

<210> SEQ ID NO 140
<211> LENGTH: 67
<212> TYPE: DNA
<213> ORGANISM: Artificial Sequence

-continued

<220> FEATURE:
<223> OTHER INFORMATION: Description of Artificial Sequence: Synthetic primer

<400> SEQUENCE: 140

ccatctcadc cctgcgtgtc tccgactcag ttattgatct ctgtacccca taaatatgta 60
caacact 67

<210> SEQ ID NO 141
<211> LENGTH: 58
<212> TYPE: DNA
<213> ORGANISM: Artificial Sequence
<220> FEATURE:
<223> OTHER INFORMATION: Description of Artificial Sequence: Synthetic primer

<400> SEQUENCE: 141

ccatctcadc cctgcgtgtc tccgactcag tacgctcgga cgccagatgg gtaagtgc 58

<210> SEQ ID NO 142
<211> LENGTH: 59
<212> TYPE: DNA
<213> ORGANISM: Artificial Sequence
<220> FEATURE:
<223> OTHER INFORMATION: Description of Artificial Sequence: Synthetic primer

<400> SEQUENCE: 142

ccatctcadc cctgcgtgtc tccgactcag caatccaagg tgcaaatcag tggtctctcc 59

<210> SEQ ID NO 143
<211> LENGTH: 58
<212> TYPE: DNA
<213> ORGANISM: Artificial Sequence
<220> FEATURE:
<223> OTHER INFORMATION: Description of Artificial Sequence: Synthetic primer

<400> SEQUENCE: 143

ccatctcadc cctgcgtgtc tccgactcag tcgtagctat ctcccttct gccacctg 58

<210> SEQ ID NO 144
<211> LENGTH: 65
<212> TYPE: DNA
<213> ORGANISM: Artificial Sequence
<220> FEATURE:
<223> OTHER INFORMATION: Description of Artificial Sequence: Synthetic primer

<400> SEQUENCE: 144

ccatctcadc cctgcgtgtc tccgactcag cgctcctgc acatgtgata cttttgggaa 60
tgaag 65

<210> SEQ ID NO 145
<211> LENGTH: 62
<212> TYPE: DNA
<213> ORGANISM: Artificial Sequence
<220> FEATURE:
<223> OTHER INFORMATION: Description of Artificial Sequence: Synthetic primer

<400> SEQUENCE: 145

ccatctcadc cctgcgtgtc tccgactcag tccgttcatt cttctgaaca ccaaattgga 60

-continued

aa 62

<210> SEQ ID NO 146
<211> LENGTH: 61
<212> TYPE: DNA
<213> ORGANISM: Artificial Sequence
<220> FEATURE:
<223> OTHER INFORMATION: Description of Artificial Sequence: Synthetic primer

<400> SEQUENCE: 146

ccatctcatc cctgcgtgtc tccgactcag cggccaggct tgtaagagc ccagagggtc 60
a 61

<210> SEQ ID NO 147
<211> LENGTH: 58
<212> TYPE: DNA
<213> ORGANISM: Artificial Sequence
<220> FEATURE:
<223> OTHER INFORMATION: Description of Artificial Sequence: Synthetic primer

<400> SEQUENCE: 147

ccatctcatc cctgcgtgtc tccgactcag caacctatct ggggcccta ctctttga 58

<210> SEQ ID NO 148
<211> LENGTH: 58
<212> TYPE: DNA
<213> ORGANISM: Artificial Sequence
<220> FEATURE:
<223> OTHER INFORMATION: Description of Artificial Sequence: Synthetic primer

<400> SEQUENCE: 148

ccatctcatc cctgcgtgtc tccgactcag cgtaatctca ctgcagccc ctacactg 58

<210> SEQ ID NO 149
<211> LENGTH: 58
<212> TYPE: DNA
<213> ORGANISM: Artificial Sequence
<220> FEATURE:
<223> OTHER INFORMATION: Description of Artificial Sequence: Synthetic primer

<400> SEQUENCE: 149

ccatctcatc cctgcgtgtc tccgactcag atatcgcgac tgactgcct ctgctgtg 58

<210> SEQ ID NO 150
<211> LENGTH: 60
<212> TYPE: DNA
<213> ORGANISM: Artificial Sequence
<220> FEATURE:
<223> OTHER INFORMATION: Description of Artificial Sequence: Synthetic primer

<400> SEQUENCE: 150

ccatctcatc cctgcgtgtc tccgactcag tcaatatctg agtgacgatg agcaggaggt 60

<210> SEQ ID NO 151
<211> LENGTH: 58
<212> TYPE: DNA
<213> ORGANISM: Artificial Sequence
<220> FEATURE:

-continued

<223> OTHER INFORMATION: Description of Artificial Sequence: Synthetic primer

<400> SEQUENCE: 151

ccatctcatc cctgcgtgtc tccgactcag atagagtata gccaatccc attgacca 58

<210> SEQ ID NO 152

<211> LENGTH: 60

<212> TYPE: DNA

<213> ORGANISM: Artificial Sequence

<220> FEATURE:

<223> OTHER INFORMATION: Description of Artificial Sequence: Synthetic primer

<400> SEQUENCE: 152

ccatctcatc cctgcgtgtc tccgactcag gcaactagtt ccaagcttcc tgaaccagac 60

<210> SEQ ID NO 153

<211> LENGTH: 58

<212> TYPE: DNA

<213> ORGANISM: Artificial Sequence

<220> FEATURE:

<223> OTHER INFORMATION: Description of Artificial Sequence: Synthetic primer

<400> SEQUENCE: 153

ccatctcatc cctgcgtgtc tccgactcag atctcgaatc ctagtggggg cattccaa 58

<210> SEQ ID NO 154

<211> LENGTH: 58

<212> TYPE: DNA

<213> ORGANISM: Artificial Sequence

<220> FEATURE:

<223> OTHER INFORMATION: Description of Artificial Sequence: Synthetic primer

<400> SEQUENCE: 154

ccatctcatc cctgcgtgtc tccgactcag ccaggagcga ctctagggcg cgtttcct 58

<210> SEQ ID NO 155

<211> LENGTH: 59

<212> TYPE: DNA

<213> ORGANISM: Artificial Sequence

<220> FEATURE:

<223> OTHER INFORMATION: Description of Artificial Sequence: Synthetic primer

<400> SEQUENCE: 155

ccatctcatc cctgcgtgtc tccgactcag atctccatcg tcagccttcc ctcgtccta 59

<210> SEQ ID NO 156

<211> LENGTH: 60

<212> TYPE: DNA

<213> ORGANISM: Artificial Sequence

<220> FEATURE:

<223> OTHER INFORMATION: Description of Artificial Sequence: Synthetic primer

<400> SEQUENCE: 156

ccatctcatc cctgcgtgtc tccgactcag ttgacgagct tccacataac tcgcttgca 60

<210> SEQ ID NO 157

<211> LENGTH: 61

<212> TYPE: DNA

-continued

<213> ORGANISM: Artificial Sequence
<220> FEATURE:
<223> OTHER INFORMATION: Description of Artificial Sequence: Synthetic primer

<400> SEQUENCE: 157

ccatctcatc cctgcgtgtc tccgactcag tactattacc gaactgtagc ccagacactg 60
c 61

<210> SEQ ID NO 158
<211> LENGTH: 63
<212> TYPE: DNA
<213> ORGANISM: Artificial Sequence
<220> FEATURE:
<223> OTHER INFORMATION: Description of Artificial Sequence: Synthetic primer

<400> SEQUENCE: 158

ccatctcatc cctgcgtgtc tccgactcag cgtcctggac acaaagctgg aaactcttcc 60
cta 63

<210> SEQ ID NO 159
<211> LENGTH: 60
<212> TYPE: DNA
<213> ORGANISM: Artificial Sequence
<220> FEATURE:
<223> OTHER INFORMATION: Description of Artificial Sequence: Synthetic primer

<400> SEQUENCE: 159

ccatctcatc cctgcgtgtc tccgactcag ctggcgctt ccaacaagcc caacaagttc 60

<210> SEQ ID NO 160
<211> LENGTH: 60
<212> TYPE: DNA
<213> ORGANISM: Artificial Sequence
<220> FEATURE:
<223> OTHER INFORMATION: Description of Artificial Sequence: Synthetic primer

<400> SEQUENCE: 160

ccatctcatc cctgcgtgtc tccgactcag gatacgtaag gaatgaccgg ctgtctgttt 60

<210> SEQ ID NO 161
<211> LENGTH: 59
<212> TYPE: DNA
<213> ORGANISM: Artificial Sequence
<220> FEATURE:
<223> OTHER INFORMATION: Description of Artificial Sequence: Synthetic primer

<400> SEQUENCE: 161

ccatctcatc cctgcgtgtc tccgactcag ctgggattaa aaagtggcac caccaacaa 59

<210> SEQ ID NO 162
<211> LENGTH: 60
<212> TYPE: DNA
<213> ORGANISM: Artificial Sequence
<220> FEATURE:
<223> OTHER INFORMATION: Description of Artificial Sequence: Synthetic primer

<400> SEQUENCE: 162

-continued

ccatctcatc cctgcgtgtc tccgactcag ttggattcgt ctgttagcgc ttcccacagt 60

<210> SEQ ID NO 163
<211> LENGTH: 61
<212> TYPE: DNA
<213> ORGANISM: Artificial Sequence
<220> FEATURE:
<223> OTHER INFORMATION: Description of Artificial Sequence: Synthetic primer

<400> SEQUENCE: 163

ccatctcatc cctgcgtgtc tccgactcag ccgtccgcta agcttctttc cacaatcctc 60

a 61

<210> SEQ ID NO 164
<211> LENGTH: 67
<212> TYPE: DNA
<213> ORGANISM: Artificial Sequence
<220> FEATURE:
<223> OTHER INFORMATION: Description of Artificial Sequence: Synthetic primer

<400> SEQUENCE: 164

ccatctcatc cctgcgtgtc tccgactcag gcgattgcaa ctgtacccca taaatatgta 60

caacact 67

<210> SEQ ID NO 165
<211> LENGTH: 58
<212> TYPE: DNA
<213> ORGANISM: Artificial Sequence
<220> FEATURE:
<223> OTHER INFORMATION: Description of Artificial Sequence: Synthetic primer

<400> SEQUENCE: 165

ccatctcatc cctgcgtgtc tccgactcag ccatgcataa cgccagatgg gtaagtgc 58

<210> SEQ ID NO 166
<211> LENGTH: 59
<212> TYPE: DNA
<213> ORGANISM: Artificial Sequence
<220> FEATURE:
<223> OTHER INFORMATION: Description of Artificial Sequence: Synthetic primer

<400> SEQUENCE: 166

ccatctcatc cctgcgtgtc tccgactcag taattgcaat tgcaaatcag tggtctctcc 59

<210> SEQ ID NO 167
<211> LENGTH: 58
<212> TYPE: DNA
<213> ORGANISM: Artificial Sequence
<220> FEATURE:
<223> OTHER INFORMATION: Description of Artificial Sequence: Synthetic primer

<400> SEQUENCE: 167

ccatctcatc cctgcgtgtc tccgactcag acgactccaa ctcccttctc gccacctg 58

<210> SEQ ID NO 168
<211> LENGTH: 65
<212> TYPE: DNA
<213> ORGANISM: Artificial Sequence

-continued

<220> FEATURE:
<223> OTHER INFORMATION: Description of Artificial Sequence: Synthetic primer
<400> SEQUENCE: 168
ccatctcatc cctgcgtgtc tccgactcag atcatgcaga acatgtgata cttttgggaa 60
tgaag 65

<210> SEQ ID NO 169
<211> LENGTH: 62
<212> TYPE: DNA
<213> ORGANISM: Artificial Sequence
<220> FEATURE:
<223> OTHER INFORMATION: Description of Artificial Sequence: Synthetic primer
<400> SEQUENCE: 169
ccatctcatc cctgcgtgtc tccgactcag aactcctaatt cttctgaaca ccaaattgga 60
aa 62

<210> SEQ ID NO 170
<211> LENGTH: 61
<212> TYPE: DNA
<213> ORGANISM: Artificial Sequence
<220> FEATURE:
<223> OTHER INFORMATION: Description of Artificial Sequence: Synthetic primer
<400> SEQUENCE: 170
ccatctcatc cctgcgtgtc tccgactcag ggatattcgt tgtaagagc ccagagggtc 60
a 61

<210> SEQ ID NO 171
<211> LENGTH: 58
<212> TYPE: DNA
<213> ORGANISM: Artificial Sequence
<220> FEATURE:
<223> OTHER INFORMATION: Description of Artificial Sequence: Synthetic primer
<400> SEQUENCE: 171
ccatctcatc cctgcgtgtc tccgactcag tcggatgact ggggcccccta ctctttga 58

<210> SEQ ID NO 172
<211> LENGTH: 58
<212> TYPE: DNA
<213> ORGANISM: Artificial Sequence
<220> FEATURE:
<223> OTHER INFORMATION: Description of Artificial Sequence: Synthetic primer
<400> SEQUENCE: 172
ccatctcatc cctgcgtgtc tccgactcag gacgcgcgag ctgcagccc ctacactg 58

<210> SEQ ID NO 173
<211> LENGTH: 58
<212> TYPE: DNA
<213> ORGANISM: Artificial Sequence
<220> FEATURE:
<223> OTHER INFORMATION: Description of Artificial Sequence: Synthetic primer
<400> SEQUENCE: 173

-continued

ccatctcattc cctgcgtgtc tccgactcag gcttagacct tgactgcct ctgctgtg 58

<210> SEQ ID NO 174
<211> LENGTH: 60
<212> TYPE: DNA
<213> ORGANISM: Artificial Sequence
<220> FEATURE:
<223> OTHER INFORMATION: Description of Artificial Sequence: Synthetic
primer

<400> SEQUENCE: 174

ccatctcattc cctgcgtgtc tccgactcag gaccaggcga agtgacgatg agcaggaggt 60

<210> SEQ ID NO 175
<211> LENGTH: 58
<212> TYPE: DNA
<213> ORGANISM: Artificial Sequence
<220> FEATURE:
<223> OTHER INFORMATION: Description of Artificial Sequence: Synthetic
primer

<400> SEQUENCE: 175

ccatctcattc cctgcgtgtc tccgactcag gctctggcgt gccaatccc attgacca 58

<210> SEQ ID NO 176
<211> LENGTH: 60
<212> TYPE: DNA
<213> ORGANISM: Artificial Sequence
<220> FEATURE:
<223> OTHER INFORMATION: Description of Artificial Sequence: Synthetic
primer

<400> SEQUENCE: 176

ccatctcattc cctgcgtgtc tccgactcag tgggccgga ccaagcttcc tgaaccagac 60

<210> SEQ ID NO 177
<211> LENGTH: 58
<212> TYPE: DNA
<213> ORGANISM: Artificial Sequence
<220> FEATURE:
<223> OTHER INFORMATION: Description of Artificial Sequence: Synthetic
primer

<400> SEQUENCE: 177

ccatctcattc cctgcgtgtc tccgactcag ctctgcgtct ctagtggggg cattccaa 58

<210> SEQ ID NO 178
<211> LENGTH: 58
<212> TYPE: DNA
<213> ORGANISM: Artificial Sequence
<220> FEATURE:
<223> OTHER INFORMATION: Description of Artificial Sequence: Synthetic
primer

<400> SEQUENCE: 178

ccatctcattc cctgcgtgtc tccgactcag ccagaagcag ctctagggcg cgtttcct 58

<210> SEQ ID NO 179
<211> LENGTH: 59
<212> TYPE: DNA
<213> ORGANISM: Artificial Sequence
<220> FEATURE:
<223> OTHER INFORMATION: Description of Artificial Sequence: Synthetic
primer

-continued

<400> SEQUENCE: 179

ccatctcatc cctgcgtgtc tccgactcag ggaaggttgc tcagccttcc ctgctcta 59

<210> SEQ ID NO 180

<211> LENGTH: 60

<212> TYPE: DNA

<213> ORGANISM: Artificial Sequence

<220> FEATURE:

<223> OTHER INFORMATION: Description of Artificial Sequence: Synthetic primer

<400> SEQUENCE: 180

ccatctcatc cctgcgtgtc tccgactcag taacggtacg tccacataac tcgcttcag 60

<210> SEQ ID NO 181

<211> LENGTH: 61

<212> TYPE: DNA

<213> ORGANISM: Artificial Sequence

<220> FEATURE:

<223> OTHER INFORMATION: Description of Artificial Sequence: Synthetic primer

<400> SEQUENCE: 181

ccatctcatc cctgcgtgtc tccgactcag ctgctcatg gaactgtagc ccagacactg 60

c 61

<210> SEQ ID NO 182

<211> LENGTH: 63

<212> TYPE: DNA

<213> ORGANISM: Artificial Sequence

<220> FEATURE:

<223> OTHER INFORMATION: Description of Artificial Sequence: Synthetic primer

<400> SEQUENCE: 182

ccatctcatc cctgcgtgtc tccgactcag actccaaggc acaaagctgg aaactcttcc 60

cta 63

<210> SEQ ID NO 183

<211> LENGTH: 60

<212> TYPE: DNA

<213> ORGANISM: Artificial Sequence

<220> FEATURE:

<223> OTHER INFORMATION: Description of Artificial Sequence: Synthetic primer

<400> SEQUENCE: 183

ccatctcatc cctgcgtgtc tccgactcag gagctgctat ccaacaagcc caacaagttc 60

<210> SEQ ID NO 184

<211> LENGTH: 60

<212> TYPE: DNA

<213> ORGANISM: Artificial Sequence

<220> FEATURE:

<223> OTHER INFORMATION: Description of Artificial Sequence: Synthetic primer

<400> SEQUENCE: 184

ccatctcatc cctgcgtgtc tccgactcag cgttgaggcc gaatgaccgg ctgtctgttt 60

<210> SEQ ID NO 185

-continued

<211> LENGTH: 59
<212> TYPE: DNA
<213> ORGANISM: Artificial Sequence
<220> FEATURE:
<223> OTHER INFORMATION: Description of Artificial Sequence: Synthetic primer

<400> SEQUENCE: 185

ccatctcatc cctgcgtgtc tccgactcag ttctggatcc aaagtggcac caccaacaa 59

<210> SEQ ID NO 186
<211> LENGTH: 60
<212> TYPE: DNA
<213> ORGANISM: Artificial Sequence
<220> FEATURE:
<223> OTHER INFORMATION: Description of Artificial Sequence: Synthetic primer

<400> SEQUENCE: 186

ccatctcatc cctgcgtgtc tccgactcag ccggattcca ctgttagcgc ttcccacagt 60

<210> SEQ ID NO 187
<211> LENGTH: 61
<212> TYPE: DNA
<213> ORGANISM: Artificial Sequence
<220> FEATURE:
<223> OTHER INFORMATION: Description of Artificial Sequence: Synthetic primer

<400> SEQUENCE: 187

ccatctcatc cctgcgtgtc tccgactcag tccatcgctt agcttctttc cacaatcctc 60
a 61

<210> SEQ ID NO 188
<211> LENGTH: 67
<212> TYPE: DNA
<213> ORGANISM: Artificial Sequence
<220> FEATURE:
<223> OTHER INFORMATION: Description of Artificial Sequence: Synthetic primer

<400> SEQUENCE: 188

ccatctcatc cctgcgtgtc tccgactcag ttactttctca ctgtacccca taaatatgta 60
caacact 67

<210> SEQ ID NO 189
<211> LENGTH: 10
<212> TYPE: DNA
<213> ORGANISM: Artificial Sequence
<220> FEATURE:
<223> OTHER INFORMATION: Description of Artificial Sequence: Synthetic oligonucleotide

<400> SEQUENCE: 189

ttaacggacg 10

<210> SEQ ID NO 190
<211> LENGTH: 10
<212> TYPE: DNA
<213> ORGANISM: Artificial Sequence
<220> FEATURE:
<223> OTHER INFORMATION: Description of Artificial Sequence: Synthetic oligonucleotide

-continued

<400> SEQUENCE: 190
tccggcttac 10

<210> SEQ ID NO 191
<211> LENGTH: 10
<212> TYPE: DNA
<213> ORGANISM: Artificial Sequence
<220> FEATURE:
<223> OTHER INFORMATION: Description of Artificial Sequence: Synthetic
oligonucleotide

<400> SEQUENCE: 191
tctcattcag 10

<210> SEQ ID NO 192
<211> LENGTH: 10
<212> TYPE: DNA
<213> ORGANISM: Artificial Sequence
<220> FEATURE:
<223> OTHER INFORMATION: Description of Artificial Sequence: Synthetic
oligonucleotide

<400> SEQUENCE: 192
gcggtcatac 10

<210> SEQ ID NO 193
<211> LENGTH: 10
<212> TYPE: DNA
<213> ORGANISM: Artificial Sequence
<220> FEATURE:
<223> OTHER INFORMATION: Description of Artificial Sequence: Synthetic
oligonucleotide

<400> SEQUENCE: 193
taggacgttc 10

<210> SEQ ID NO 194
<211> LENGTH: 10
<212> TYPE: DNA
<213> ORGANISM: Artificial Sequence
<220> FEATURE:
<223> OTHER INFORMATION: Description of Artificial Sequence: Synthetic
oligonucleotide

<400> SEQUENCE: 194
acgacgcaac 10

<210> SEQ ID NO 195
<211> LENGTH: 10
<212> TYPE: DNA
<213> ORGANISM: Artificial Sequence
<220> FEATURE:
<223> OTHER INFORMATION: Description of Artificial Sequence: Synthetic
oligonucleotide

<400> SEQUENCE: 195
cttctcggac 10

<210> SEQ ID NO 196
<211> LENGTH: 10
<212> TYPE: DNA
<213> ORGANISM: Artificial Sequence
<220> FEATURE:
<223> OTHER INFORMATION: Description of Artificial Sequence: Synthetic

-continued

oligonucleotide

<400> SEQUENCE: 196

cattgccgtt 10

<210> SEQ ID NO 197
<211> LENGTH: 10
<212> TYPE: DNA
<213> ORGANISM: Artificial Sequence
<220> FEATURE:
<223> OTHER INFORMATION: Description of Artificial Sequence: Synthetic
oligonucleotide

<400> SEQUENCE: 197

cgagccagaa 10

<210> SEQ ID NO 198
<211> LENGTH: 10
<212> TYPE: DNA
<213> ORGANISM: Artificial Sequence
<220> FEATURE:
<223> OTHER INFORMATION: Description of Artificial Sequence: Synthetic
oligonucleotide

<400> SEQUENCE: 198

tgaggacggc 10

<210> SEQ ID NO 199
<211> LENGTH: 10
<212> TYPE: DNA
<213> ORGANISM: Artificial Sequence
<220> FEATURE:
<223> OTHER INFORMATION: Description of Artificial Sequence: Synthetic
oligonucleotide

<400> SEQUENCE: 199

gcctgcgag 10

<210> SEQ ID NO 200
<211> LENGTH: 10
<212> TYPE: DNA
<213> ORGANISM: Artificial Sequence
<220> FEATURE:
<223> OTHER INFORMATION: Description of Artificial Sequence: Synthetic
oligonucleotide

<400> SEQUENCE: 200

gttgacgtct 10

<210> SEQ ID NO 201
<211> LENGTH: 10
<212> TYPE: DNA
<213> ORGANISM: Artificial Sequence
<220> FEATURE:
<223> OTHER INFORMATION: Description of Artificial Sequence: Synthetic
oligonucleotide

<400> SEQUENCE: 201

gagatcgatt 10

<210> SEQ ID NO 202
<211> LENGTH: 10
<212> TYPE: DNA
<213> ORGANISM: Artificial Sequence

-continued

<220> FEATURE:
<223> OTHER INFORMATION: Description of Artificial Sequence: Synthetic
oligonucleotide

<400> SEQUENCE: 202

agttcgagcc 10

<210> SEQ ID NO 203
<211> LENGTH: 10
<212> TYPE: DNA
<213> ORGANISM: Artificial Sequence
<220> FEATURE:
<223> OTHER INFORMATION: Description of Artificial Sequence: Synthetic
oligonucleotide

<400> SEQUENCE: 203

ctcagggtca 10

<210> SEQ ID NO 204
<211> LENGTH: 10
<212> TYPE: DNA
<213> ORGANISM: Artificial Sequence
<220> FEATURE:
<223> OTHER INFORMATION: Description of Artificial Sequence: Synthetic
oligonucleotide

<400> SEQUENCE: 204

ggcaatataa 10

<210> SEQ ID NO 205
<211> LENGTH: 10
<212> TYPE: DNA
<213> ORGANISM: Artificial Sequence
<220> FEATURE:
<223> OTHER INFORMATION: Description of Artificial Sequence: Synthetic
oligonucleotide

<400> SEQUENCE: 205

ggtactcatg 10

<210> SEQ ID NO 206
<211> LENGTH: 10
<212> TYPE: DNA
<213> ORGANISM: Artificial Sequence
<220> FEATURE:
<223> OTHER INFORMATION: Description of Artificial Sequence: Synthetic
oligonucleotide

<400> SEQUENCE: 206

tctggttcaa 10

<210> SEQ ID NO 207
<211> LENGTH: 10
<212> TYPE: DNA
<213> ORGANISM: Artificial Sequence
<220> FEATURE:
<223> OTHER INFORMATION: Description of Artificial Sequence: Synthetic
oligonucleotide

<400> SEQUENCE: 207

gtcctataag 10

<210> SEQ ID NO 208
<211> LENGTH: 10

-continued

<212> TYPE: DNA
<213> ORGANISM: Artificial Sequence
<220> FEATURE:
<223> OTHER INFORMATION: Description of Artificial Sequence: Synthetic
oligonucleotide

<400> SEQUENCE: 208

gtcagcctcc 10

<210> SEQ ID NO 209
<211> LENGTH: 10
<212> TYPE: DNA
<213> ORGANISM: Artificial Sequence
<220> FEATURE:
<223> OTHER INFORMATION: Description of Artificial Sequence: Synthetic
oligonucleotide

<400> SEQUENCE: 209

ttcaagctcg 10

<210> SEQ ID NO 210
<211> LENGTH: 10
<212> TYPE: DNA
<213> ORGANISM: Artificial Sequence
<220> FEATURE:
<223> OTHER INFORMATION: Description of Artificial Sequence: Synthetic
oligonucleotide

<400> SEQUENCE: 210

gtaccagcgc 10

<210> SEQ ID NO 211
<211> LENGTH: 10
<212> TYPE: DNA
<213> ORGANISM: Artificial Sequence
<220> FEATURE:
<223> OTHER INFORMATION: Description of Artificial Sequence: Synthetic
oligonucleotide

<400> SEQUENCE: 211

tcctattcgg 10

<210> SEQ ID NO 212
<211> LENGTH: 10
<212> TYPE: DNA
<213> ORGANISM: Artificial Sequence
<220> FEATURE:
<223> OTHER INFORMATION: Description of Artificial Sequence: Synthetic
oligonucleotide

<400> SEQUENCE: 212

gccagcgatt 10

<210> SEQ ID NO 213
<211> LENGTH: 10
<212> TYPE: DNA
<213> ORGANISM: Artificial Sequence
<220> FEATURE:
<223> OTHER INFORMATION: Description of Artificial Sequence: Synthetic
oligonucleotide

<400> SEQUENCE: 213

acctagactg 10

-continued

<210> SEQ ID NO 214
<211> LENGTH: 10
<212> TYPE: DNA
<213> ORGANISM: Artificial Sequence
<220> FEATURE:
<223> OTHER INFORMATION: Description of Artificial Sequence: Synthetic
oligonucleotide

<400> SEQUENCE: 214

actggttcgc 10

<210> SEQ ID NO 215
<211> LENGTH: 10
<212> TYPE: DNA
<213> ORGANISM: Artificial Sequence
<220> FEATURE:
<223> OTHER INFORMATION: Description of Artificial Sequence: Synthetic
oligonucleotide

<400> SEQUENCE: 215

ccatattagg 10

<210> SEQ ID NO 216
<211> LENGTH: 10
<212> TYPE: DNA
<213> ORGANISM: Artificial Sequence
<220> FEATURE:
<223> OTHER INFORMATION: Description of Artificial Sequence: Synthetic
oligonucleotide

<400> SEQUENCE: 216

gctcgtcagc 10

<210> SEQ ID NO 217
<211> LENGTH: 10
<212> TYPE: DNA
<213> ORGANISM: Artificial Sequence
<220> FEATURE:
<223> OTHER INFORMATION: Description of Artificial Sequence: Synthetic
oligonucleotide

<400> SEQUENCE: 217

cgtaatgacg 10

<210> SEQ ID NO 218
<211> LENGTH: 10
<212> TYPE: DNA
<213> ORGANISM: Artificial Sequence
<220> FEATURE:
<223> OTHER INFORMATION: Description of Artificial Sequence: Synthetic
oligonucleotide

<400> SEQUENCE: 218

ccggcgctga 10

<210> SEQ ID NO 219
<211> LENGTH: 10
<212> TYPE: DNA
<213> ORGANISM: Artificial Sequence
<220> FEATURE:
<223> OTHER INFORMATION: Description of Artificial Sequence: Synthetic
oligonucleotide

<400> SEQUENCE: 219

cgcgaagata 10

-continued

<210> SEQ ID NO 220
<211> LENGTH: 10
<212> TYPE: DNA
<213> ORGANISM: Artificial Sequence
<220> FEATURE:
<223> OTHER INFORMATION: Description of Artificial Sequence: Synthetic
oligonucleotide

<400> SEQUENCE: 220

gaaccgcaga 10

<210> SEQ ID NO 221
<211> LENGTH: 10
<212> TYPE: DNA
<213> ORGANISM: Artificial Sequence
<220> FEATURE:
<223> OTHER INFORMATION: Description of Artificial Sequence: Synthetic
oligonucleotide

<400> SEQUENCE: 221

ttggcagaga 10

<210> SEQ ID NO 222
<211> LENGTH: 10
<212> TYPE: DNA
<213> ORGANISM: Artificial Sequence
<220> FEATURE:
<223> OTHER INFORMATION: Description of Artificial Sequence: Synthetic
oligonucleotide

<400> SEQUENCE: 222

gcattctctgc 10

<210> SEQ ID NO 223
<211> LENGTH: 10
<212> TYPE: DNA
<213> ORGANISM: Artificial Sequence
<220> FEATURE:
<223> OTHER INFORMATION: Description of Artificial Sequence: Synthetic
oligonucleotide

<400> SEQUENCE: 223

ttggaccgca 10

<210> SEQ ID NO 224
<211> LENGTH: 10
<212> TYPE: DNA
<213> ORGANISM: Artificial Sequence
<220> FEATURE:
<223> OTHER INFORMATION: Description of Artificial Sequence: Synthetic
oligonucleotide

<400> SEQUENCE: 224

gcagaacgtc 10

<210> SEQ ID NO 225
<211> LENGTH: 10
<212> TYPE: DNA
<213> ORGANISM: Artificial Sequence
<220> FEATURE:
<223> OTHER INFORMATION: Description of Artificial Sequence: Synthetic
oligonucleotide

<400> SEQUENCE: 225

-continued

aacttcgagc 10

<210> SEQ ID NO 226
<211> LENGTH: 10
<212> TYPE: DNA
<213> ORGANISM: Artificial Sequence
<220> FEATURE:
<223> OTHER INFORMATION: Description of Artificial Sequence: Synthetic
oligonucleotide

<400> SEQUENCE: 226

gctcctagag 10

<210> SEQ ID NO 227
<211> LENGTH: 10
<212> TYPE: DNA
<213> ORGANISM: Artificial Sequence
<220> FEATURE:
<223> OTHER INFORMATION: Description of Artificial Sequence: Synthetic
oligonucleotide

<400> SEQUENCE: 227

tatctagctt 10

<210> SEQ ID NO 228
<211> LENGTH: 10
<212> TYPE: DNA
<213> ORGANISM: Artificial Sequence
<220> FEATURE:
<223> OTHER INFORMATION: Description of Artificial Sequence: Synthetic
oligonucleotide

<400> SEQUENCE: 228

gagtattggc 10

<210> SEQ ID NO 229
<211> LENGTH: 10
<212> TYPE: DNA
<213> ORGANISM: Artificial Sequence
<220> FEATURE:
<223> OTHER INFORMATION: Description of Artificial Sequence: Synthetic
oligonucleotide

<400> SEQUENCE: 229

cctgagctca 10

<210> SEQ ID NO 230
<211> LENGTH: 10
<212> TYPE: DNA
<213> ORGANISM: Artificial Sequence
<220> FEATURE:
<223> OTHER INFORMATION: Description of Artificial Sequence: Synthetic
oligonucleotide

<400> SEQUENCE: 230

caggcgagta 10

<210> SEQ ID NO 231
<211> LENGTH: 10
<212> TYPE: DNA
<213> ORGANISM: Artificial Sequence
<220> FEATURE:
<223> OTHER INFORMATION: Description of Artificial Sequence: Synthetic
oligonucleotide

-continued

<400> SEQUENCE: 231

gcaggcagag 10

<210> SEQ ID NO 232

<211> LENGTH: 10

<212> TYPE: DNA

<213> ORGANISM: Artificial Sequence

<220> FEATURE:

<223> OTHER INFORMATION: Description of Artificial Sequence: Synthetic
oligonucleotide

<400> SEQUENCE: 232

gcgtcgatac 10

<210> SEQ ID NO 233

<211> LENGTH: 10

<212> TYPE: DNA

<213> ORGANISM: Artificial Sequence

<220> FEATURE:

<223> OTHER INFORMATION: Description of Artificial Sequence: Synthetic
oligonucleotide

<400> SEQUENCE: 233

cgatgattat 10

<210> SEQ ID NO 234

<211> LENGTH: 10

<212> TYPE: DNA

<213> ORGANISM: Artificial Sequence

<220> FEATURE:

<223> OTHER INFORMATION: Description of Artificial Sequence: Synthetic
oligonucleotide

<400> SEQUENCE: 234

gacggctggc 10

<210> SEQ ID NO 235

<211> LENGTH: 10

<212> TYPE: DNA

<213> ORGANISM: Artificial Sequence

<220> FEATURE:

<223> OTHER INFORMATION: Description of Artificial Sequence: Synthetic
oligonucleotide

<400> SEQUENCE: 235

ggagcctgag 10

<210> SEQ ID NO 236

<211> LENGTH: 10

<212> TYPE: DNA

<213> ORGANISM: Artificial Sequence

<220> FEATURE:

<223> OTHER INFORMATION: Description of Artificial Sequence: Synthetic
oligonucleotide

<400> SEQUENCE: 236

cctgactgct 10

<210> SEQ ID NO 237

<211> LENGTH: 10

<212> TYPE: DNA

<213> ORGANISM: Artificial Sequence

<220> FEATURE:

-continued

<223> OTHER INFORMATION: Description of Artificial Sequence: Synthetic
oligonucleotide

<400> SEQUENCE: 237

acggctgacg 10

<210> SEQ ID NO 238
<211> LENGTH: 10
<212> TYPE: DNA
<213> ORGANISM: Artificial Sequence
<220> FEATURE:
<223> OTHER INFORMATION: Description of Artificial Sequence: Synthetic
oligonucleotide

<400> SEQUENCE: 238

taaccatagc 10

<210> SEQ ID NO 239
<211> LENGTH: 10
<212> TYPE: DNA
<213> ORGANISM: Artificial Sequence
<220> FEATURE:
<223> OTHER INFORMATION: Description of Artificial Sequence: Synthetic
oligonucleotide

<400> SEQUENCE: 239

tcttgcttc 10

<210> SEQ ID NO 240
<211> LENGTH: 10
<212> TYPE: DNA
<213> ORGANISM: Artificial Sequence
<220> FEATURE:
<223> OTHER INFORMATION: Description of Artificial Sequence: Synthetic
oligonucleotide

<400> SEQUENCE: 240

ttcttagatt 10

<210> SEQ ID NO 241
<211> LENGTH: 10
<212> TYPE: DNA
<213> ORGANISM: Artificial Sequence
<220> FEATURE:
<223> OTHER INFORMATION: Description of Artificial Sequence: Synthetic
oligonucleotide

<400> SEQUENCE: 241

tcatctcatt 10

<210> SEQ ID NO 242
<211> LENGTH: 10
<212> TYPE: DNA
<213> ORGANISM: Artificial Sequence
<220> FEATURE:
<223> OTHER INFORMATION: Description of Artificial Sequence: Synthetic
oligonucleotide

<400> SEQUENCE: 242

tctccgctcg 10

<210> SEQ ID NO 243
<211> LENGTH: 10
<212> TYPE: DNA

-continued

<213> ORGANISM: Artificial Sequence
<220> FEATURE:
<223> OTHER INFORMATION: Description of Artificial Sequence: Synthetic
oligonucleotide

<400> SEQUENCE: 243

tgccatatgc 10

<210> SEQ ID NO 244
<211> LENGTH: 10
<212> TYPE: DNA
<213> ORGANISM: Artificial Sequence
<220> FEATURE:
<223> OTHER INFORMATION: Description of Artificial Sequence: Synthetic
oligonucleotide

<400> SEQUENCE: 244

taaggcctct 10

<210> SEQ ID NO 245
<211> LENGTH: 10
<212> TYPE: DNA
<213> ORGANISM: Artificial Sequence
<220> FEATURE:
<223> OTHER INFORMATION: Description of Artificial Sequence: Synthetic
oligonucleotide

<400> SEQUENCE: 245

gagtaggccg 10

<210> SEQ ID NO 246
<211> LENGTH: 10
<212> TYPE: DNA
<213> ORGANISM: Artificial Sequence
<220> FEATURE:
<223> OTHER INFORMATION: Description of Artificial Sequence: Synthetic
oligonucleotide

<400> SEQUENCE: 246

gcaataagct 10

<210> SEQ ID NO 247
<211> LENGTH: 10
<212> TYPE: DNA
<213> ORGANISM: Artificial Sequence
<220> FEATURE:
<223> OTHER INFORMATION: Description of Artificial Sequence: Synthetic
oligonucleotide

<400> SEQUENCE: 247

ggcgttgcaa 10

<210> SEQ ID NO 248
<211> LENGTH: 10
<212> TYPE: DNA
<213> ORGANISM: Artificial Sequence
<220> FEATURE:
<223> OTHER INFORMATION: Description of Artificial Sequence: Synthetic
oligonucleotide

<400> SEQUENCE: 248

ccaagaagcg 10

<210> SEQ ID NO 249

-continued

<211> LENGTH: 10
<212> TYPE: DNA
<213> ORGANISM: Artificial Sequence
<220> FEATURE:
<223> OTHER INFORMATION: Description of Artificial Sequence: Synthetic
oligonucleotide

<400> SEQUENCE: 249

ggttacctcg 10

<210> SEQ ID NO 250
<211> LENGTH: 10
<212> TYPE: DNA
<213> ORGANISM: Artificial Sequence
<220> FEATURE:
<223> OTHER INFORMATION: Description of Artificial Sequence: Synthetic
oligonucleotide

<400> SEQUENCE: 250

ctccgcctta 10

<210> SEQ ID NO 251
<211> LENGTH: 10
<212> TYPE: DNA
<213> ORGANISM: Artificial Sequence
<220> FEATURE:
<223> OTHER INFORMATION: Description of Artificial Sequence: Synthetic
oligonucleotide

<400> SEQUENCE: 251

ctccagagat 10

<210> SEQ ID NO 252
<211> LENGTH: 10
<212> TYPE: DNA
<213> ORGANISM: Artificial Sequence
<220> FEATURE:
<223> OTHER INFORMATION: Description of Artificial Sequence: Synthetic
oligonucleotide

<400> SEQUENCE: 252

gtcgaggtag 10

<210> SEQ ID NO 253
<211> LENGTH: 10
<212> TYPE: DNA
<213> ORGANISM: Artificial Sequence
<220> FEATURE:
<223> OTHER INFORMATION: Description of Artificial Sequence: Synthetic
oligonucleotide

<400> SEQUENCE: 253

tatggacctg 10

<210> SEQ ID NO 254
<211> LENGTH: 10
<212> TYPE: DNA
<213> ORGANISM: Artificial Sequence
<220> FEATURE:
<223> OTHER INFORMATION: Description of Artificial Sequence: Synthetic
oligonucleotide

<400> SEQUENCE: 254

tacctgctag 10

-continued

<210> SEQ ID NO 255
<211> LENGTH: 10
<212> TYPE: DNA
<213> ORGANISM: Artificial Sequence
<220> FEATURE:
<223> OTHER INFORMATION: Description of Artificial Sequence: Synthetic
oligonucleotide

<400> SEQUENCE: 255

ccgcgaccga 10

<210> SEQ ID NO 256
<211> LENGTH: 10
<212> TYPE: DNA
<213> ORGANISM: Artificial Sequence
<220> FEATURE:
<223> OTHER INFORMATION: Description of Artificial Sequence: Synthetic
oligonucleotide

<400> SEQUENCE: 256

gttgaacgtt 10

<210> SEQ ID NO 257
<211> LENGTH: 10
<212> TYPE: DNA
<213> ORGANISM: Artificial Sequence
<220> FEATURE:
<223> OTHER INFORMATION: Description of Artificial Sequence: Synthetic
oligonucleotide

<400> SEQUENCE: 257

tgccaacgca 10

<210> SEQ ID NO 258
<211> LENGTH: 10
<212> TYPE: DNA
<213> ORGANISM: Artificial Sequence
<220> FEATURE:
<223> OTHER INFORMATION: Description of Artificial Sequence: Synthetic
oligonucleotide

<400> SEQUENCE: 258

ggattgacct 10

<210> SEQ ID NO 259
<211> LENGTH: 10
<212> TYPE: DNA
<213> ORGANISM: Artificial Sequence
<220> FEATURE:
<223> OTHER INFORMATION: Description of Artificial Sequence: Synthetic
oligonucleotide

<400> SEQUENCE: 259

ggacggattc 10

<210> SEQ ID NO 260
<211> LENGTH: 10
<212> TYPE: DNA
<213> ORGANISM: Artificial Sequence
<220> FEATURE:
<223> OTHER INFORMATION: Description of Artificial Sequence: Synthetic
oligonucleotide

<400> SEQUENCE: 260

-continued

tcctccgtcg	10
 <210> SEQ ID NO 261 <211> LENGTH: 10 <212> TYPE: DNA <213> ORGANISM: Artificial Sequence <220> FEATURE: <223> OTHER INFORMATION: Description of Artificial Sequence: Synthetic oligonucleotide <400> SEQUENCE: 261	
agttcatggt	10
 <210> SEQ ID NO 262 <211> LENGTH: 10 <212> TYPE: DNA <213> ORGANISM: Artificial Sequence <220> FEATURE: <223> OTHER INFORMATION: Description of Artificial Sequence: Synthetic oligonucleotide <400> SEQUENCE: 262	
tatccattcc	10
 <210> SEQ ID NO 263 <211> LENGTH: 10 <212> TYPE: DNA <213> ORGANISM: Artificial Sequence <220> FEATURE: <223> OTHER INFORMATION: Description of Artificial Sequence: Synthetic oligonucleotide <400> SEQUENCE: 263	
ggagagcgcg	10
 <210> SEQ ID NO 264 <211> LENGTH: 10 <212> TYPE: DNA <213> ORGANISM: Artificial Sequence <220> FEATURE: <223> OTHER INFORMATION: Description of Artificial Sequence: Synthetic oligonucleotide <400> SEQUENCE: 264	
cggaccttgg	10
 <210> SEQ ID NO 265 <211> LENGTH: 10 <212> TYPE: DNA <213> ORGANISM: Artificial Sequence <220> FEATURE: <223> OTHER INFORMATION: Description of Artificial Sequence: Synthetic oligonucleotide <400> SEQUENCE: 265	
ggcaatctcc	10
 <210> SEQ ID NO 266 <211> LENGTH: 10 <212> TYPE: DNA <213> ORGANISM: Artificial Sequence <220> FEATURE: <223> OTHER INFORMATION: Description of Artificial Sequence: Synthetic oligonucleotide	

-continued

<400> SEQUENCE: 266

aggattgatt 10

<210> SEQ ID NO 267

<211> LENGTH: 10

<212> TYPE: DNA

<213> ORGANISM: Artificial Sequence

<220> FEATURE:

<223> OTHER INFORMATION: Description of Artificial Sequence: Synthetic
oligonucleotide

<400> SEQUENCE: 267

gccgttgctt 10

<210> SEQ ID NO 268

<211> LENGTH: 10

<212> TYPE: DNA

<213> ORGANISM: Artificial Sequence

<220> FEATURE:

<223> OTHER INFORMATION: Description of Artificial Sequence: Synthetic
oligonucleotide

<400> SEQUENCE: 268

aagtacgtcg 10

<210> SEQ ID NO 269

<211> LENGTH: 10

<212> TYPE: DNA

<213> ORGANISM: Artificial Sequence

<220> FEATURE:

<223> OTHER INFORMATION: Description of Artificial Sequence: Synthetic
oligonucleotide

<400> SEQUENCE: 269

tggcttaagg 10

<210> SEQ ID NO 270

<211> LENGTH: 10

<212> TYPE: DNA

<213> ORGANISM: Artificial Sequence

<220> FEATURE:

<223> OTHER INFORMATION: Description of Artificial Sequence: Synthetic
oligonucleotide

<400> SEQUENCE: 270

ctcttcaga 10

<210> SEQ ID NO 271

<211> LENGTH: 10

<212> TYPE: DNA

<213> ORGANISM: Artificial Sequence

<220> FEATURE:

<223> OTHER INFORMATION: Description of Artificial Sequence: Synthetic
oligonucleotide

<400> SEQUENCE: 271

cgttcttcaa 10

<210> SEQ ID NO 272

<211> LENGTH: 10

<212> TYPE: DNA

<213> ORGANISM: Artificial Sequence

<220> FEATURE:

<223> OTHER INFORMATION: Description of Artificial Sequence: Synthetic

-continued

oligonucleotide

<400> SEQUENCE: 272

caacggctgc 10

<210> SEQ ID NO 273
<211> LENGTH: 10
<212> TYPE: DNA
<213> ORGANISM: Artificial Sequence
<220> FEATURE:
<223> OTHER INFORMATION: Description of Artificial Sequence: Synthetic
oligonucleotide

<400> SEQUENCE: 273

gcaagtaacc 10

<210> SEQ ID NO 274
<211> LENGTH: 10
<212> TYPE: DNA
<213> ORGANISM: Artificial Sequence
<220> FEATURE:
<223> OTHER INFORMATION: Description of Artificial Sequence: Synthetic
oligonucleotide

<400> SEQUENCE: 274

gttcatagtc 10

<210> SEQ ID NO 275
<211> LENGTH: 10
<212> TYPE: DNA
<213> ORGANISM: Artificial Sequence
<220> FEATURE:
<223> OTHER INFORMATION: Description of Artificial Sequence: Synthetic
oligonucleotide

<400> SEQUENCE: 275

acggcgagcc 10

<210> SEQ ID NO 276
<211> LENGTH: 10
<212> TYPE: DNA
<213> ORGANISM: Artificial Sequence
<220> FEATURE:
<223> OTHER INFORMATION: Description of Artificial Sequence: Synthetic
oligonucleotide

<400> SEQUENCE: 276

gtatggtcgg 10

<210> SEQ ID NO 277
<211> LENGTH: 10
<212> TYPE: DNA
<213> ORGANISM: Artificial Sequence
<220> FEATURE:
<223> OTHER INFORMATION: Description of Artificial Sequence: Synthetic
oligonucleotide

<400> SEQUENCE: 277

tcggttatcc 10

<210> SEQ ID NO 278
<211> LENGTH: 10
<212> TYPE: DNA
<213> ORGANISM: Artificial Sequence

-continued

<220> FEATURE:
<223> OTHER INFORMATION: Description of Artificial Sequence: Synthetic
oligonucleotide

<400> SEQUENCE: 278

gcggtcgata 10

<210> SEQ ID NO 279
<211> LENGTH: 10
<212> TYPE: DNA
<213> ORGANISM: Artificial Sequence
<220> FEATURE:
<223> OTHER INFORMATION: Description of Artificial Sequence: Synthetic
oligonucleotide

<400> SEQUENCE: 279

tcctcagtat 10

<210> SEQ ID NO 280
<211> LENGTH: 10
<212> TYPE: DNA
<213> ORGANISM: Artificial Sequence
<220> FEATURE:
<223> OTHER INFORMATION: Description of Artificial Sequence: Synthetic
oligonucleotide

<400> SEQUENCE: 280

accgttcctg 10

<210> SEQ ID NO 281
<211> LENGTH: 10
<212> TYPE: DNA
<213> ORGANISM: Artificial Sequence
<220> FEATURE:
<223> OTHER INFORMATION: Description of Artificial Sequence: Synthetic
oligonucleotide

<400> SEQUENCE: 281

gcctgctctt 10

<210> SEQ ID NO 282
<211> LENGTH: 10
<212> TYPE: DNA
<213> ORGANISM: Artificial Sequence
<220> FEATURE:
<223> OTHER INFORMATION: Description of Artificial Sequence: Synthetic
oligonucleotide

<400> SEQUENCE: 282

agcgtaacca 10

<210> SEQ ID NO 283
<211> LENGTH: 10
<212> TYPE: DNA
<213> ORGANISM: Artificial Sequence
<220> FEATURE:
<223> OTHER INFORMATION: Description of Artificial Sequence: Synthetic
oligonucleotide

<400> SEQUENCE: 283

ttgcctgatg 10

<210> SEQ ID NO 284
<211> LENGTH: 10

-continued

<212> TYPE: DNA
<213> ORGANISM: Artificial Sequence
<220> FEATURE:
<223> OTHER INFORMATION: Description of Artificial Sequence: Synthetic
oligonucleotide

<400> SEQUENCE: 284

ttattgatct 10

<210> SEQ ID NO 285
<211> LENGTH: 10
<212> TYPE: DNA
<213> ORGANISM: Artificial Sequence
<220> FEATURE:
<223> OTHER INFORMATION: Description of Artificial Sequence: Synthetic
oligonucleotide

<400> SEQUENCE: 285

tacgctcgga 10

<210> SEQ ID NO 286
<211> LENGTH: 10
<212> TYPE: DNA
<213> ORGANISM: Artificial Sequence
<220> FEATURE:
<223> OTHER INFORMATION: Description of Artificial Sequence: Synthetic
oligonucleotide

<400> SEQUENCE: 286

caatccaagg 10

<210> SEQ ID NO 287
<211> LENGTH: 10
<212> TYPE: DNA
<213> ORGANISM: Artificial Sequence
<220> FEATURE:
<223> OTHER INFORMATION: Description of Artificial Sequence: Synthetic
oligonucleotide

<400> SEQUENCE: 287

tcgtagctat 10

<210> SEQ ID NO 288
<211> LENGTH: 10
<212> TYPE: DNA
<213> ORGANISM: Artificial Sequence
<220> FEATURE:
<223> OTHER INFORMATION: Description of Artificial Sequence: Synthetic
oligonucleotide

<400> SEQUENCE: 288

cgctcatcg 10

<210> SEQ ID NO 289
<211> LENGTH: 10
<212> TYPE: DNA
<213> ORGANISM: Artificial Sequence
<220> FEATURE:
<223> OTHER INFORMATION: Description of Artificial Sequence: Synthetic
oligonucleotide

<400> SEQUENCE: 289

tccgttcatt 10

-continued

<210> SEQ ID NO 290
<211> LENGTH: 10
<212> TYPE: DNA
<213> ORGANISM: Artificial Sequence
<220> FEATURE:
<223> OTHER INFORMATION: Description of Artificial Sequence: Synthetic
oligonucleotide

<400> SEQUENCE: 290

cggccaggct 10

<210> SEQ ID NO 291
<211> LENGTH: 10
<212> TYPE: DNA
<213> ORGANISM: Artificial Sequence
<220> FEATURE:
<223> OTHER INFORMATION: Description of Artificial Sequence: Synthetic
oligonucleotide

<400> SEQUENCE: 291

caacctatct 10

<210> SEQ ID NO 292
<211> LENGTH: 10
<212> TYPE: DNA
<213> ORGANISM: Artificial Sequence
<220> FEATURE:
<223> OTHER INFORMATION: Description of Artificial Sequence: Synthetic
oligonucleotide

<400> SEQUENCE: 292

cgtaatctca 10

<210> SEQ ID NO 293
<211> LENGTH: 10
<212> TYPE: DNA
<213> ORGANISM: Artificial Sequence
<220> FEATURE:
<223> OTHER INFORMATION: Description of Artificial Sequence: Synthetic
oligonucleotide

<400> SEQUENCE: 293

atatcgcgac 10

<210> SEQ ID NO 294
<211> LENGTH: 10
<212> TYPE: DNA
<213> ORGANISM: Artificial Sequence
<220> FEATURE:
<223> OTHER INFORMATION: Description of Artificial Sequence: Synthetic
oligonucleotide

<400> SEQUENCE: 294

tcaatatctg 10

<210> SEQ ID NO 295
<211> LENGTH: 10
<212> TYPE: DNA
<213> ORGANISM: Artificial Sequence
<220> FEATURE:
<223> OTHER INFORMATION: Description of Artificial Sequence: Synthetic
oligonucleotide

<400> SEQUENCE: 295

atagagtata 10

-continued

<210> SEQ ID NO 296
<211> LENGTH: 10
<212> TYPE: DNA
<213> ORGANISM: Artificial Sequence
<220> FEATURE:
<223> OTHER INFORMATION: Description of Artificial Sequence: Synthetic
oligonucleotide

<400> SEQUENCE: 296

gcaactagtt 10

<210> SEQ ID NO 297
<211> LENGTH: 10
<212> TYPE: DNA
<213> ORGANISM: Artificial Sequence
<220> FEATURE:
<223> OTHER INFORMATION: Description of Artificial Sequence: Synthetic
oligonucleotide

<400> SEQUENCE: 297

atctcgaatc 10

<210> SEQ ID NO 298
<211> LENGTH: 10
<212> TYPE: DNA
<213> ORGANISM: Artificial Sequence
<220> FEATURE:
<223> OTHER INFORMATION: Description of Artificial Sequence: Synthetic
oligonucleotide

<400> SEQUENCE: 298

ccaggagcga 10

<210> SEQ ID NO 299
<211> LENGTH: 10
<212> TYPE: DNA
<213> ORGANISM: Artificial Sequence
<220> FEATURE:
<223> OTHER INFORMATION: Description of Artificial Sequence: Synthetic
oligonucleotide

<400> SEQUENCE: 299

atctccatcg 10

<210> SEQ ID NO 300
<211> LENGTH: 10
<212> TYPE: DNA
<213> ORGANISM: Artificial Sequence
<220> FEATURE:
<223> OTHER INFORMATION: Description of Artificial Sequence: Synthetic
oligonucleotide

<400> SEQUENCE: 300

ttgacgagct 10

<210> SEQ ID NO 301
<211> LENGTH: 10
<212> TYPE: DNA
<213> ORGANISM: Artificial Sequence
<220> FEATURE:
<223> OTHER INFORMATION: Description of Artificial Sequence: Synthetic
oligonucleotide

<400> SEQUENCE: 301

-continued

tactattacc 10

<210> SEQ ID NO 302
<211> LENGTH: 10
<212> TYPE: DNA
<213> ORGANISM: Artificial Sequence
<220> FEATURE:
<223> OTHER INFORMATION: Description of Artificial Sequence: Synthetic
oligonucleotide

<400> SEQUENCE: 302

cgtcctggac 10

<210> SEQ ID NO 303
<211> LENGTH: 10
<212> TYPE: DNA
<213> ORGANISM: Artificial Sequence
<220> FEATURE:
<223> OTHER INFORMATION: Description of Artificial Sequence: Synthetic
oligonucleotide

<400> SEQUENCE: 303

ctcggcgctt 10

<210> SEQ ID NO 304
<211> LENGTH: 10
<212> TYPE: DNA
<213> ORGANISM: Artificial Sequence
<220> FEATURE:
<223> OTHER INFORMATION: Description of Artificial Sequence: Synthetic
oligonucleotide

<400> SEQUENCE: 304

gatacgtaag 10

<210> SEQ ID NO 305
<211> LENGTH: 10
<212> TYPE: DNA
<213> ORGANISM: Artificial Sequence
<220> FEATURE:
<223> OTHER INFORMATION: Description of Artificial Sequence: Synthetic
oligonucleotide

<400> SEQUENCE: 305

ctcggattaa 10

<210> SEQ ID NO 306
<211> LENGTH: 10
<212> TYPE: DNA
<213> ORGANISM: Artificial Sequence
<220> FEATURE:
<223> OTHER INFORMATION: Description of Artificial Sequence: Synthetic
oligonucleotide

<400> SEQUENCE: 306

ttggattcgt 10

<210> SEQ ID NO 307
<211> LENGTH: 10
<212> TYPE: DNA
<213> ORGANISM: Artificial Sequence
<220> FEATURE:
<223> OTHER INFORMATION: Description of Artificial Sequence: Synthetic
oligonucleotide

-continued

<400> SEQUENCE: 307

ccgtccgcta 10

<210> SEQ ID NO 308

<211> LENGTH: 10

<212> TYPE: DNA

<213> ORGANISM: Artificial Sequence

<220> FEATURE:

<223> OTHER INFORMATION: Description of Artificial Sequence: Synthetic
oligonucleotide

<400> SEQUENCE: 308

gcgattgcaa 10

<210> SEQ ID NO 309

<211> LENGTH: 10

<212> TYPE: DNA

<213> ORGANISM: Artificial Sequence

<220> FEATURE:

<223> OTHER INFORMATION: Description of Artificial Sequence: Synthetic
oligonucleotide

<400> SEQUENCE: 309

ccatgcataa 10

<210> SEQ ID NO 310

<211> LENGTH: 10

<212> TYPE: DNA

<213> ORGANISM: Artificial Sequence

<220> FEATURE:

<223> OTHER INFORMATION: Description of Artificial Sequence: Synthetic
oligonucleotide

<400> SEQUENCE: 310

taattgcaat 10

<210> SEQ ID NO 311

<211> LENGTH: 10

<212> TYPE: DNA

<213> ORGANISM: Artificial Sequence

<220> FEATURE:

<223> OTHER INFORMATION: Description of Artificial Sequence: Synthetic
oligonucleotide

<400> SEQUENCE: 311

acgactccaa 10

<210> SEQ ID NO 312

<211> LENGTH: 10

<212> TYPE: DNA

<213> ORGANISM: Artificial Sequence

<220> FEATURE:

<223> OTHER INFORMATION: Description of Artificial Sequence: Synthetic
oligonucleotide

<400> SEQUENCE: 312

atcatgcaga 10

<210> SEQ ID NO 313

<211> LENGTH: 10

<212> TYPE: DNA

<213> ORGANISM: Artificial Sequence

<220> FEATURE:

-continued

<223> OTHER INFORMATION: Description of Artificial Sequence: Synthetic oligonucleotide

<400> SEQUENCE: 313

aactcctaatt 10

<210> SEQ ID NO 314
<211> LENGTH: 10
<212> TYPE: DNA
<213> ORGANISM: Artificial Sequence
<220> FEATURE:
<223> OTHER INFORMATION: Description of Artificial Sequence: Synthetic oligonucleotide

<400> SEQUENCE: 314

ggatattcgt 10

<210> SEQ ID NO 315
<211> LENGTH: 10
<212> TYPE: DNA
<213> ORGANISM: Artificial Sequence
<220> FEATURE:
<223> OTHER INFORMATION: Description of Artificial Sequence: Synthetic oligonucleotide

<400> SEQUENCE: 315

tcggatgact 10

<210> SEQ ID NO 316
<211> LENGTH: 10
<212> TYPE: DNA
<213> ORGANISM: Artificial Sequence
<220> FEATURE:
<223> OTHER INFORMATION: Description of Artificial Sequence: Synthetic oligonucleotide

<400> SEQUENCE: 316

gacgcgcgag 10

<210> SEQ ID NO 317
<211> LENGTH: 10
<212> TYPE: DNA
<213> ORGANISM: Artificial Sequence
<220> FEATURE:
<223> OTHER INFORMATION: Description of Artificial Sequence: Synthetic oligonucleotide

<400> SEQUENCE: 317

gcctagacct 10

<210> SEQ ID NO 318
<211> LENGTH: 10
<212> TYPE: DNA
<213> ORGANISM: Artificial Sequence
<220> FEATURE:
<223> OTHER INFORMATION: Description of Artificial Sequence: Synthetic oligonucleotide

<400> SEQUENCE: 318

gaccaggcga 10

<210> SEQ ID NO 319
<211> LENGTH: 10
<212> TYPE: DNA

-continued

<213> ORGANISM: Artificial Sequence
<220> FEATURE:
<223> OTHER INFORMATION: Description of Artificial Sequence: Synthetic
oligonucleotide

<400> SEQUENCE: 319

gctctggcgt 10

<210> SEQ ID NO 320
<211> LENGTH: 10
<212> TYPE: DNA
<213> ORGANISM: Artificial Sequence
<220> FEATURE:
<223> OTHER INFORMATION: Description of Artificial Sequence: Synthetic
oligonucleotide

<400> SEQUENCE: 320

tggtccgga 10

<210> SEQ ID NO 321
<211> LENGTH: 10
<212> TYPE: DNA
<213> ORGANISM: Artificial Sequence
<220> FEATURE:
<223> OTHER INFORMATION: Description of Artificial Sequence: Synthetic
oligonucleotide

<400> SEQUENCE: 321

ctctgcgtct 10

<210> SEQ ID NO 322
<211> LENGTH: 10
<212> TYPE: DNA
<213> ORGANISM: Artificial Sequence
<220> FEATURE:
<223> OTHER INFORMATION: Description of Artificial Sequence: Synthetic
oligonucleotide

<400> SEQUENCE: 322

ccagaagcag 10

<210> SEQ ID NO 323
<211> LENGTH: 10
<212> TYPE: DNA
<213> ORGANISM: Artificial Sequence
<220> FEATURE:
<223> OTHER INFORMATION: Description of Artificial Sequence: Synthetic
oligonucleotide

<400> SEQUENCE: 323

ggaaggttgc 10

<210> SEQ ID NO 324
<211> LENGTH: 10
<212> TYPE: DNA
<213> ORGANISM: Artificial Sequence
<220> FEATURE:
<223> OTHER INFORMATION: Description of Artificial Sequence: Synthetic
oligonucleotide

<400> SEQUENCE: 324

taacggtacg 10

<210> SEQ ID NO 325

-continued

<211> LENGTH: 10
<212> TYPE: DNA
<213> ORGANISM: Artificial Sequence
<220> FEATURE:
<223> OTHER INFORMATION: Description of Artificial Sequence: Synthetic
oligonucleotide

<400> SEQUENCE: 325

ctcgctcatg 10

<210> SEQ ID NO 326
<211> LENGTH: 10
<212> TYPE: DNA
<213> ORGANISM: Artificial Sequence
<220> FEATURE:
<223> OTHER INFORMATION: Description of Artificial Sequence: Synthetic
oligonucleotide

<400> SEQUENCE: 326

actccaaggc 10

<210> SEQ ID NO 327
<211> LENGTH: 10
<212> TYPE: DNA
<213> ORGANISM: Artificial Sequence
<220> FEATURE:
<223> OTHER INFORMATION: Description of Artificial Sequence: Synthetic
oligonucleotide

<400> SEQUENCE: 327

gagctgctat 10

<210> SEQ ID NO 328
<211> LENGTH: 10
<212> TYPE: DNA
<213> ORGANISM: Artificial Sequence
<220> FEATURE:
<223> OTHER INFORMATION: Description of Artificial Sequence: Synthetic
oligonucleotide

<400> SEQUENCE: 328

cgttgaggcc 10

<210> SEQ ID NO 329
<211> LENGTH: 10
<212> TYPE: DNA
<213> ORGANISM: Artificial Sequence
<220> FEATURE:
<223> OTHER INFORMATION: Description of Artificial Sequence: Synthetic
oligonucleotide

<400> SEQUENCE: 329

ttctggatcc 10

<210> SEQ ID NO 330
<211> LENGTH: 10
<212> TYPE: DNA
<213> ORGANISM: Artificial Sequence
<220> FEATURE:
<223> OTHER INFORMATION: Description of Artificial Sequence: Synthetic
oligonucleotide

<400> SEQUENCE: 330

ccggattcca 10

-continued

<210> SEQ ID NO 331
<211> LENGTH: 10
<212> TYPE: DNA
<213> ORGANISM: Artificial Sequence
<220> FEATURE:
<223> OTHER INFORMATION: Description of Artificial Sequence: Synthetic
oligonucleotide

<400> SEQUENCE: 331

tccatcgctt 10

<210> SEQ ID NO 332
<211> LENGTH: 10
<212> TYPE: DNA
<213> ORGANISM: Artificial Sequence
<220> FEATURE:
<223> OTHER INFORMATION: Description of Artificial Sequence: Synthetic
oligonucleotide

<400> SEQUENCE: 332

ttactttctca 10

<210> SEQ ID NO 333
<211> LENGTH: 73
<212> TYPE: DNA
<213> ORGANISM: Homo sapiens

<400> SEQUENCE: 333

aaatgcactc tgatatgggc ctaatcatat aatgccaatc tcccagggtc caaagggtga 60

tgccggtgca cct 73

<210> SEQ ID NO 334
<211> LENGTH: 9
<212> TYPE: PRT
<213> ORGANISM: Homo sapiens

<400> SEQUENCE: 334

Gly Pro Lys Gly Asp Ala Gly Ala Pro
1 5

<210> SEQ ID NO 335
<211> LENGTH: 15
<212> TYPE: DNA
<213> ORGANISM: Homo sapiens

<400> SEQUENCE: 335

tttggcagtg tcttc 15

<210> SEQ ID NO 336
<211> LENGTH: 100
<212> TYPE: DNA
<213> ORGANISM: Homo sapiens

<400> SEQUENCE: 336

ggcaactctg tagatctagg atgttcggtc agcgattcat tcttttcttc taaagatttt 60

tccaatttgg tgttcagaag attatttttt agtaagttat 100

<210> SEQ ID NO 337
<211> LENGTH: 34
<212> TYPE: PRT
<213> ORGANISM: Homo sapiens

-continued

<400> SEQUENCE: 337

```

Asn Asn Leu Leu Lys Asn Asn Leu Leu Asn Thr Lys Leu Glu Lys Ser
1      5      10      15
Leu Glu Glu Lys Asn Glu Ser Leu Thr Glu His Pro Arg Ser Thr Glu
      20      25      30
Leu Pro

```

<210> SEQ ID NO 338

<211> LENGTH: 36

<212> TYPE: DNA

<213> ORGANISM: Homo sapiens

<400> SEQUENCE: 338

ttgcgccccg ggetgccgcg cctgcgggt gcccc

36

<210> SEQ ID NO 339

<211> LENGTH: 12

<212> TYPE: PRT

<213> ORGANISM: Homo sapiens

<400> SEQUENCE: 339

```

Gly Gly Ser Arg Gln Gly Gly Gln Pro Gly Ala Gln
1      5      10

```

1. A method for determining alternate allele frequency, the method comprising:

- a) performing two or more parallel amplification reactions on a single sample, thereby generating overlapping amplicons, wherein each amplification reaction comprises a unique pair of forward and reverse primers, wherein the forward or reverse primer comprises an index sequence, and wherein the forward and reverse primers comprise different adapter sequences;
- b) sequencing the overlapping amplicons to produce sequence reads;
- c) segregating the sequencing reads into bins by index sequence; and
- d) detecting the presence or absence of one or more genetic variants within sequencing reads within a bin, wherein the frequency of detection of the variant determines the alternate allele frequency.

2. A method for determining alternate allele frequency, the method comprising:

- a) performing three amplification reactions on a single sample, thereby generating three overlapping amplicons, wherein each amplification reaction comprises a unique pair of forward and reverse primers, wherein each primer comprises a nucleic acid sequence complementary to a portion of a target nucleic acid sequence, wherein the forward or reverse primer comprises an index sequence, and wherein the forward and reverse primers comprise different adapter sequences at or near the 5' terminus of the primer and upstream of the sequence complementary to the target, and wherein at least one adapter sequence is complementary to a nucleic acid sequence used in sequencing;
- b) sequencing the overlapping amplicons to produce sequence reads;

c) segregating the sequencing reads into bins by index sequence; and

- d) detecting the presence or absence of one or more genetic variants within sequencing reads within a bin, wherein the frequency of detection of the variant determines the alternate allele frequency.

3. A method for determining alternate allele frequency, the method comprising:

- a) performing three amplification reactions on a single sample, thereby generating three overlapping amplicons, wherein each amplification reaction comprises a unique pair of forward and reverse primers, wherein the forward or reverse primer comprises an index sequence and/or a unique molecular identifier (UMI); and each primer comprises
 - i. a nucleotide sequence complementary to a portion of a target nucleic acid sequence;
 - ii. an adapter at or near its 5' terminus, wherein the adapter is upstream of the sequence complementary to the target and wherein the forward and reverse primers comprise different adapter sequences, wherein at least one adapter sequence is complementary to a nucleic acid sequence used in sequencing;
- b) sequencing the overlapping amplicons to produce sequence reads;
- c) segregating the sequencing reads into bins by index sequence;
- d) detecting the UMI and removing duplicate reads from the bin, wherein the detecting can be simultaneous with step c or subsequent to step c; and
- e) detecting the presence or absence of one or more genetic variants within sequencing reads within a bin, wherein the frequency of detection of the variant determines the alternate allele frequency.

4. The method of claim 1 further comprising pooling the amplicons prior to sequencing.

5. The method of claim 1, wherein sequencing the amplicons comprises contacting the amplicons with a nucleic acid complementary to the adapter sequence.

6. The method of claim 1, wherein the amplicons comprise a nucleotide having a label, optionally wherein the label is biotin.

7. (canceled)

8. The method of claim 6 further comprising contacting the label with a capture agent that specifically binds the label.

9. The method of claim 1 further comprising enzymatically digesting the primers.

10. The method of claim 1 further comprising amplifying the amplicons, thereby generating enriched populations of amplicons.

11. The method of claim 1, wherein the genetic variation to be detected is known or unknown.

12. The method of claim 1, wherein the genetic variant has an alternate allele fraction of at least 0.1%.

13. The method of claim 1, wherein the genetic variant has an alternate allele fraction of at least 0.025%.

14. The method of claim 1, wherein the genetic variant is a mosaic variant.

15. The method of claim 1, wherein detection of the genetic variant identifies the presence of a disease or a predisposition to a disease in a subject from whom the sample was derived.

16. The method of claim 15, wherein the disease is cancer.

17. The method claim 1, wherein the sample comprises circulating tumor cells or cell free DNA.

18. The method of claim 1, wherein the genetic variant originated from a somatic event or a germline event.

19. The method of claim 15, wherein the alternate allele frequency is compared to the allele frequency of a reference sample to determine if the subject's disease is progressing, regressing, or in remission.

20. The method of claim 1 further comprising averaging the alternate allele frequencies determined for each bin.

21. The method of claim 20 further comprising determining the error rate of the nucleic acid sequences flanking the alternate allele.

* * * * *