

Give Feedback

i-CAT FLX V-Series

Cone Beam 3D Imaging.The i-CAT FLX V-Series offers three maximum fields of view (FOV): V8, V10, and V17. Within these sizes, you have the option to scale or collimate the scan height to capture only the area of interest per the patient's immediate need. The system also provides enhanced low-dose imaging options, dedicated traditional 2D panoramic capabilities, and the feature-rich treatment planning software.

Osseo 2022
U.S. Patent No. 6,381,301

Want to see i-CAT FLX V-Series in action?

Schedule a Demo

It's hard to choose when all three models have class-leading performance.

i-CAT FLX V8 – With 8 cm x 5 cm to 8 cm x 8 cm FOVs, this system accommodates those dentists who want to capture select treatment areas for survey, implant, oral surgery, and endodontics.

i-CAT FLX V10 – The V10 offers a maximum and scalable FOV of up to 16 cm D x 10 cm H, helpful when you expand into larger cases and airway or TMJ analysis.

i-CAT FLX V17 – With a maximum FOV of up to 23 cm D x 17 cm H, the V17 is used most often by orthodontics, oral maxillofacial surgeons, and oral radiologists.

Give Feedback

Proactive Servicing through DEXIS Connect Pro


Keep your digital equipment operating at 100%

Learn more



i-CAT Highlights

Give Feedback

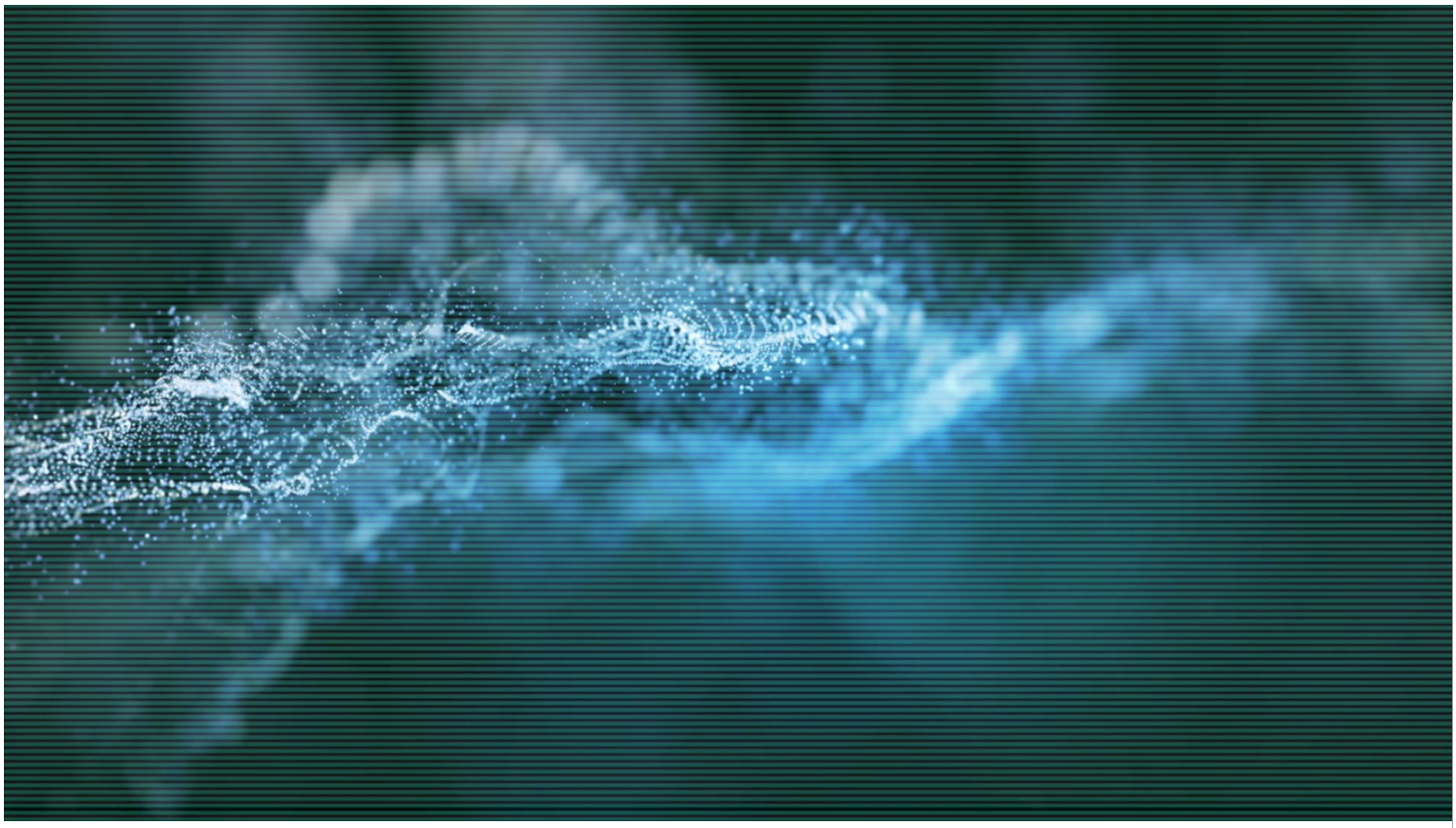

60 DAYS

SIXTY DAY

SATISFACTION

GUARANTEED

Give Feedback

At DEXIS™, we want you and your team to be confident and comfortable with your new technology. With our 60 Day Satisfaction Guarantee, we stand by you, and we stand behind the quality of our products.

See the Guarantee



With the purchase of any new CBCT equipment, DEXIS and BeamReaders are teaming up to offer 30 radiology reports in your first 90 days of ownership. Whether you are new to CBCT or an imaging veteran, enhance your radiological and diagnostic capabilities as you incorporate your new CBCT into your practice.

Get the Full Report

Give Feedback

Overview



Give Feedback

i-CAT FLX V17

Designed to accommodate the widest range of clinical applications, the V17 is the most flexible of the V-Series options. The V17 has a maximum and

scalable FOV of up to 23 cm x 17 cm. The V17 is used most often by orthodontists, oral maxillofacial surgeons, and oral radiologists; however, this system is appropriate for any dentist who wishes to view the entire oral-facial complex.

Applications
Applications for the i-CAT FLX V17 offer the most flexibility and image sizes for all 3D needs including orthodontics, orthognathics, airway analysis, and oral surgery.

Compatible with Tx STUDIO
Tx STUDIO, powered by Anatomage, is an integral part of the fast i-CAT workflow and provides the power of multiple software systems combined into one simple-to-use solution.

i-CAT FLX V10

The i-CAT FLX V10 combines highly precise cone beam 3D technology with flexible planning and treatment tools. i-CAT offers a full suite of solutions t meet your practice's needs. You can care for your patients with greater confidence and control than ever before.

The i-CAT FLX V10 has all the functionality of the V8, such as clinical control, powerful comprehensive treatment tools, and more advanced procedures with greater predictability. The difference between the V8 and V10 primarily relates to the available fields of view. V10 scan sizes range from an adjustable 8cm x 5cm to 16cm x 10cm.

Applications
Applications for the i-CAT FLX V10 include implant, most OMS, periodontics, upper/lower airway analysis, TMJ, and pediatrics.

Give Feedback

Compatible with Tx STUDIO

Tx STUDIO, powered by Anatomage, is an integral part of the fast i-CAT workflow and provides the power of multiple software systems combined into one simple-to-use solution.



Give Feedback



Give Feedback

i-CAT FLX V8

The V8 represents your gateway into the world of industry-leading 3D CBCT imaging at an affordable price. With 8cm x 5cm and 8cm x 8cm FOVs, clinicians are able to capture specific areas of the anatomy instead of the entire oral-maxillofacial complex accomplished in one scan. Should the need arise for more coverage of the anatomy, upgrading to the V10 is

simple. All V-Series systems use the same capture process, and the scan is viewed within the same 3D software.

Applications

Applications for the i-CAT FLX V8 include single implant, general practice, and endodontics.

Compatible with Tx STUDIO

Tx STUDIO, powered by Anatomage, is an integral part of the fast i-CAT workflow and provides the power of multiple software systems combined into one simple-to-use solution.

Features



Give Feedback

Low-Dose Options

i-CAT QuickScan and QuickScan+ protocols allow you to take complete 3D images at a radiation dose comparable to a 2D panoramic image*Utilizing the i-CAT FLX QuickScan+ exposure tool protocol. Use of lower dosage imaging may only be suitable for certain diagnostic tasks. Image quality is proportional to dose. i-CAT FLX offers a variety of exposure protocols allowing clinicians to adjust dosage to specific diagnostic needs.




Give Feedback

Seated Position

i-CAT FLX is one of the few machines that provides a seated position for scanning. Ergonomic Stability System (ESS) allows for easy, seated patient

positioning, designed to minimize patient movement, and avoid unnecessary retakes and additional radiation.

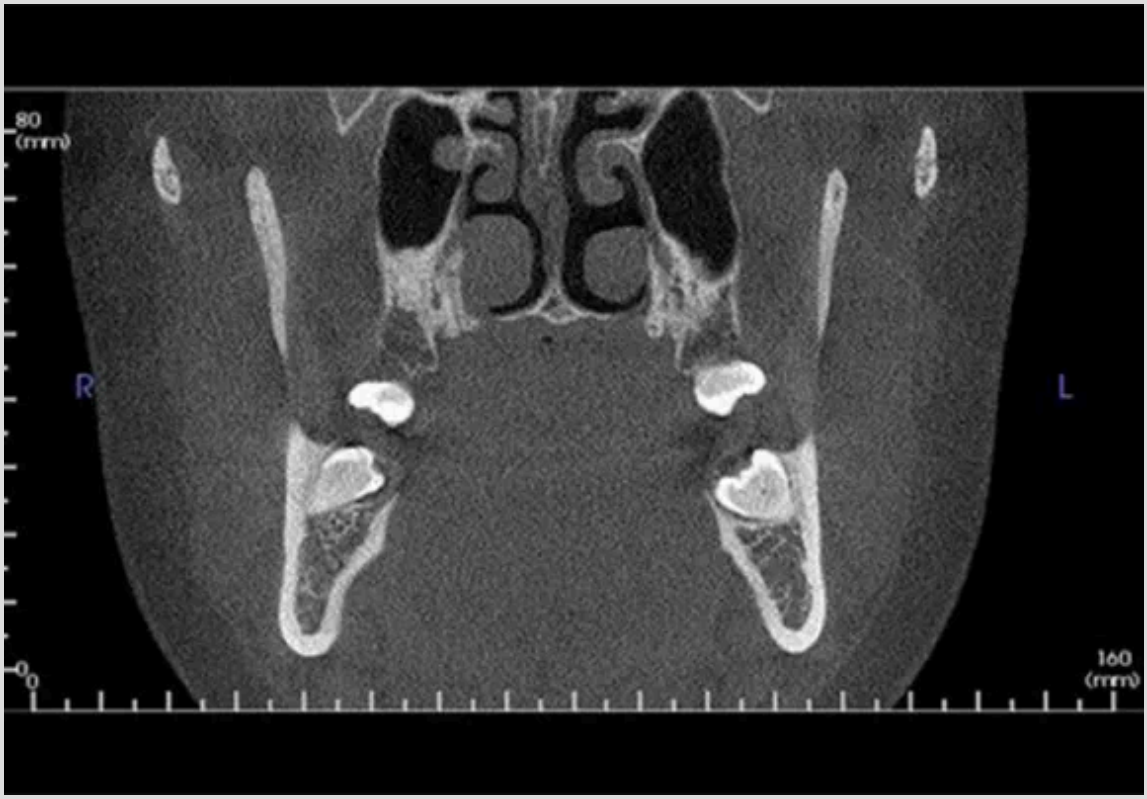

Give Feedback

Open Platform

i-CAT FLX allows for an open platform when using your DICOM files. These files are not ‘locked’, giving you the freedom to use i-CAT-captured files in a wide array of programs including CAD/CAM applications and other dental equipment.



Traditional 2D Panoramics

Your “2-in-1” i-CAT FLX V-Series takes a traditional 2D pan using the same high-quality sensor that is used to acquire 3D scans. Using your system to capture BOTH 2D and 3D imaging assists you with your ALARA goals and helps you transition to 3D imaging at your own pace.


Give Feedback



Full 24.2 cm x 19.3 cm Sensor

All i-CAT FLX V-Series models come equipped with a full 24.2 cm x 19.3 cm sensor capable of yielding a full ceph height 3D scan (13 cm x 16 cm) in just 4.8 seconds. This allows owners the flexible freedom to upgrade their field-of-view without the need to replace the unit itself. The solution is simply an in-office software upgrade that opens the field-of-view to the next V-Series configuration.

Give Feedback

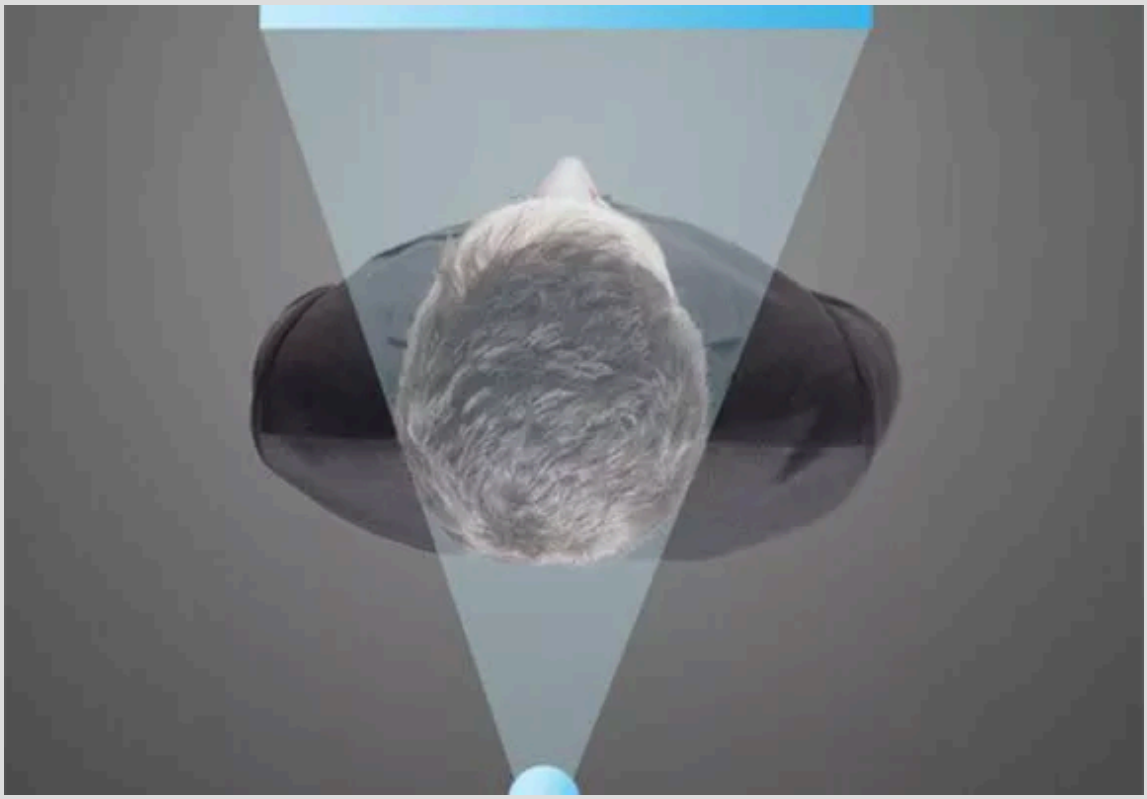

No Offset Scanning

With our PureScan™ technology and full-size sensor, i-CAT FLX V-Series scanners capture a 3D scan without compromising image quality. This full-beam scanning maximizes the use of the large sensor to provide more anatomical information at a higher quality compared to other smaller sensors.

Unlike the i-CAT FLX V-Series, some 3D CBCT imaging machines use a smaller sensor, which only captures half of the volume at a time requiring images to be stitched together. Offset or half-beam scanning allows a smaller sensor to shift and scan a slightly larger volume than its size,

Give Feedback

although the cost associated is the sacrifice of image quality and

Clinical Indications

Your dental practice is unique—that's why you need a flexible 3D solution that works with your practice and provides the control you need to deliver optimum care to your patients.

The i-CAT FLX V-series offers effective treatment planning and predictability to help general dentists and specialists alike by providing fast sizable scans, low-dose options, and dedicated specialty tools within Tx STUDIO, powered by Anatomage.

Give Feedback

Implantology


Using the i-CAT FLX V8, V10 or V17, clinicians can place and restore with accuracy and confidence.

Treat patients with greater surgical predictability and confident outcomes using i-CAT's 3D treatment planning tools. Use i-CAT's high resolution, volumetric images and complete 3D views for a more thorough analysis of bone structure and tooth orientation. Collect precise data and map an

entire course of treatment for surgical placement of the implant and abutment, all the way to final restoration.






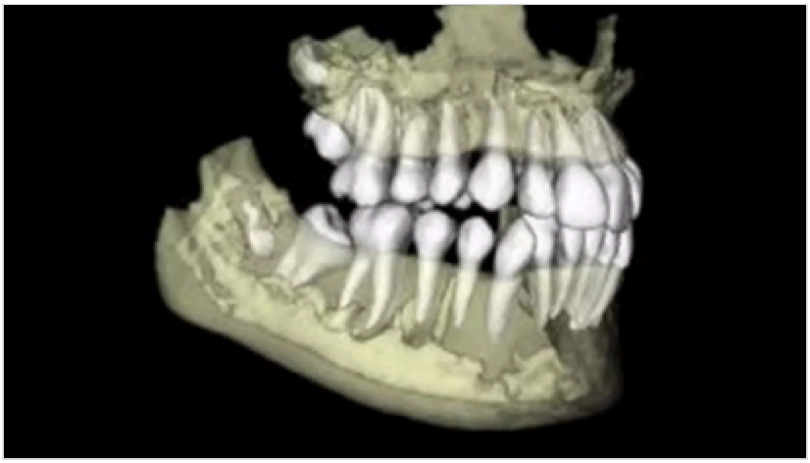


Give Feedback



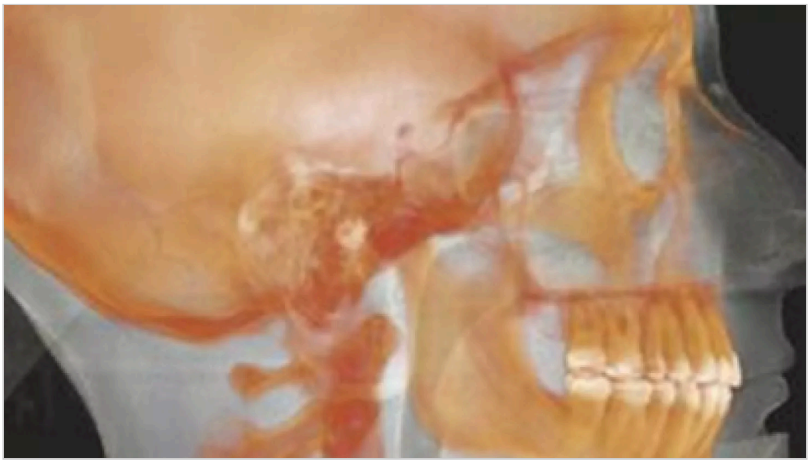




Give Feedback

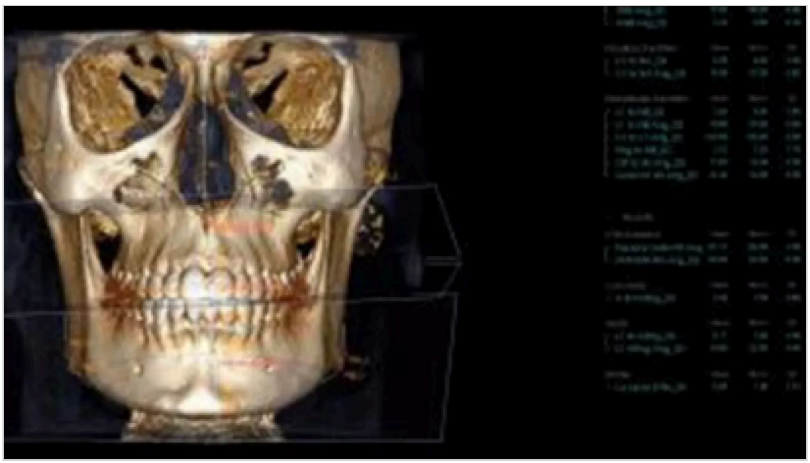



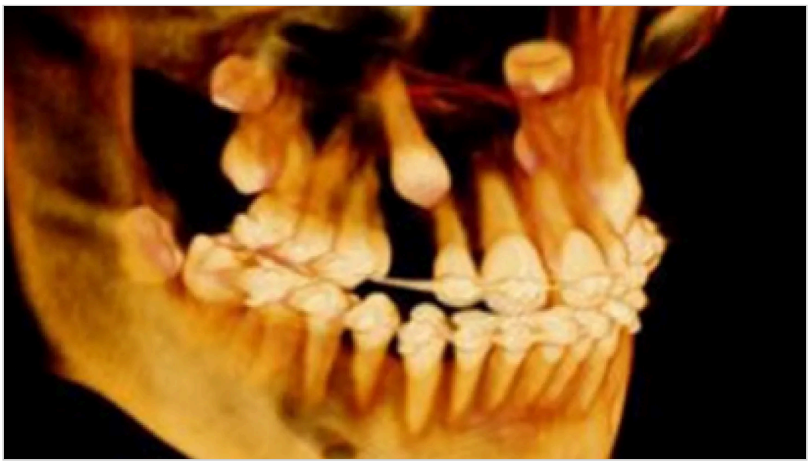

Give Feedback

Oral and Maxillofacial Surgery

Tx STUDIO treatment planning software, powered by Anatomage, can assist in identifying deformities, such as cysts, tumors, lesions, and changes of the


jaw, to avoid potential surgical complications. Determine precise position of impacted teeth within the alveolar bone, as well as their proximity to adjacent teeth and vital structures, such as the nerve canal, sinus walls, and cortical borders.

Using the i-CAT FLX V8, V10 or V17, clinicians can place and restore with accuracy and confidence.

Treat patients with greater surgical predictability and confident outcomes using i-CAT's 3D treatment planning tools.

Use i-CAT's high resolution, volumetric images and complete 3D views for a more thorough analysis of bone structure and tooth orientation.

Collect precise data and map an entire course of treatment for surgical placement of the implant and abutment, all the way to final restoration.

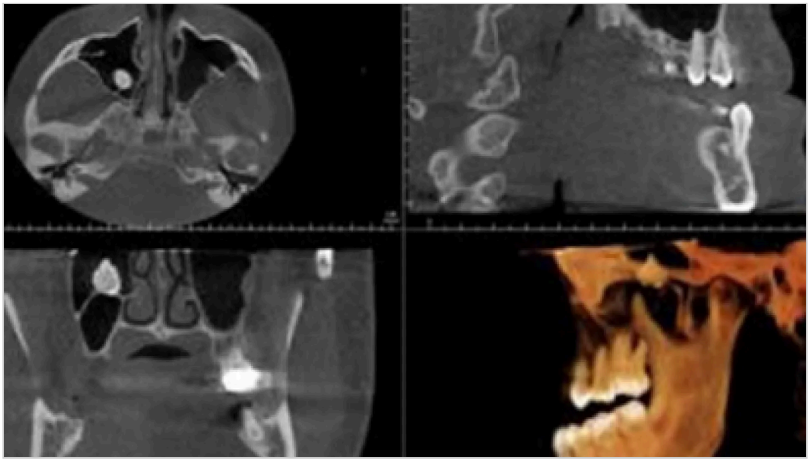

Give Feedback










Give Feedback




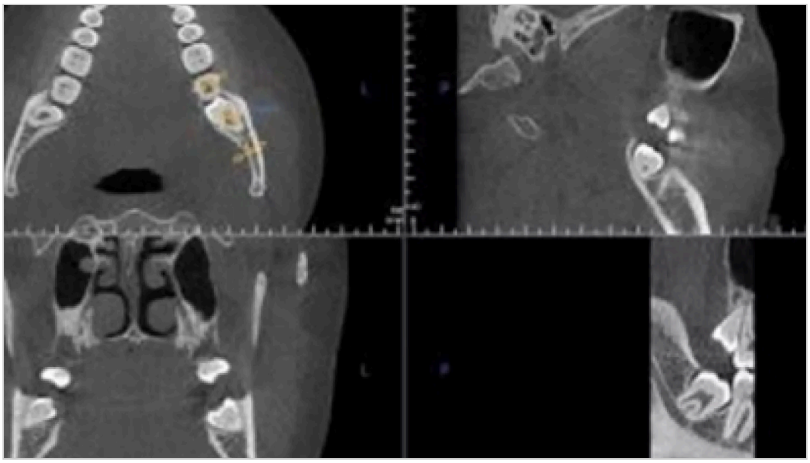

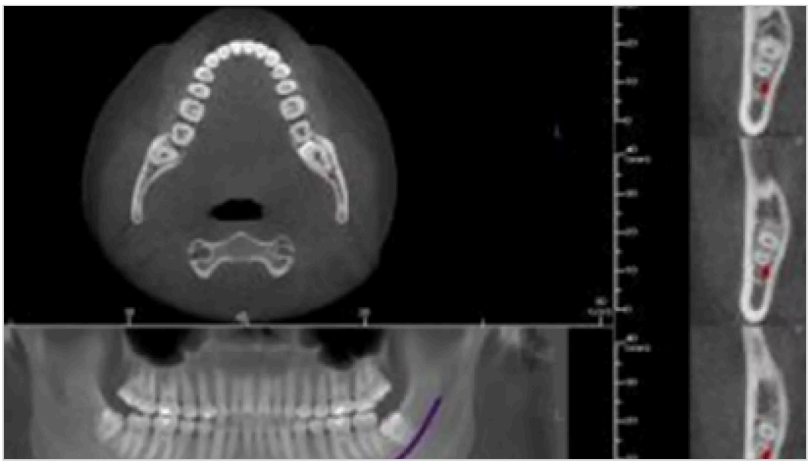

Give Feedback

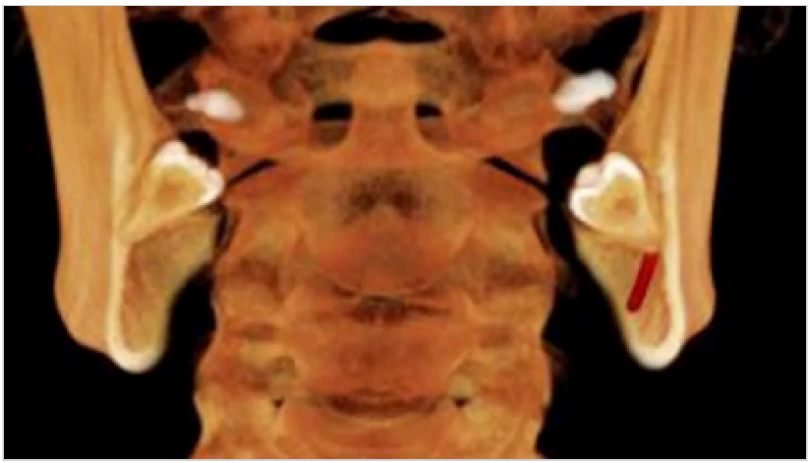

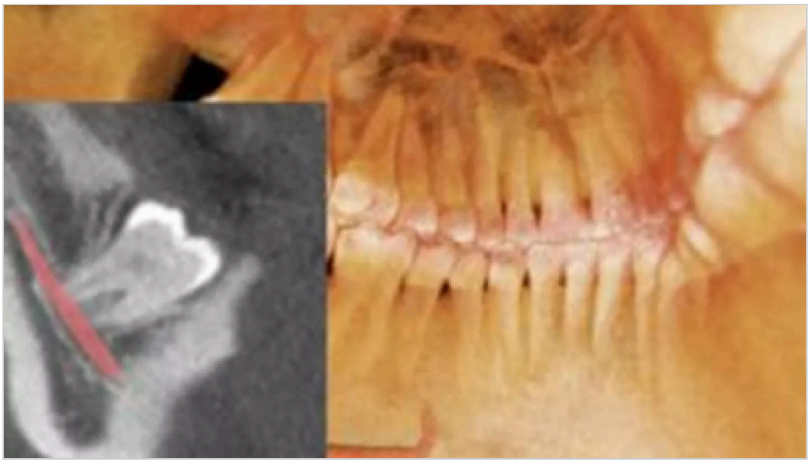

Give Feedback

Orthodontics

With the I-CAT FLX V17, clinicians can understand exact tooth position and the relationship of anatomy so they can map the most effective, least invasive treatment plan for the best alignment. Correct root angulations and find supernumerary teeth and their exact locations to enhance communication with oral surgeons. Additional treatment modules in the 3D treatment planning software, including 3D cephalometric analysis, virtual studies and impressionless models, make planning even more powerful.

Enhanced communication with oral surgeons
Clinicians can utilize the image data provided by i-CAT FLX V17 and the 3D treatment planning software to correct root angulations and find exact

locations of supernumerary teeth. Exploratory surgery may be prevented based upon the clinician's better understanding of the patient's conditions.

Additional Treatment Modules

Planning is even more powerful with 3D cephalometric analysis, virtual studies, and impression-less models. Incorporate CBCT scans in digital orthodontic treatment plans and personalized appliances.






Give Feedback






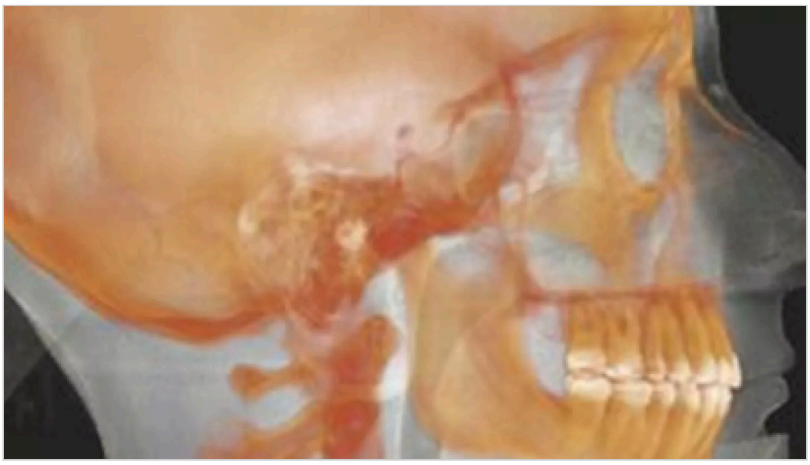


Give Feedback



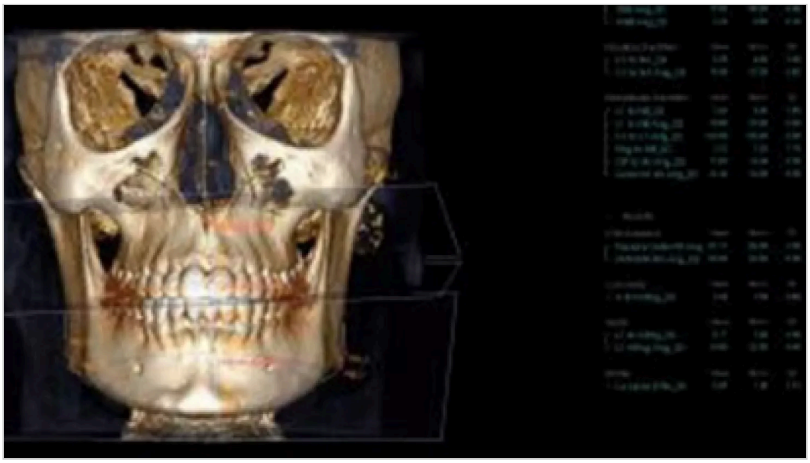



Give Feedback

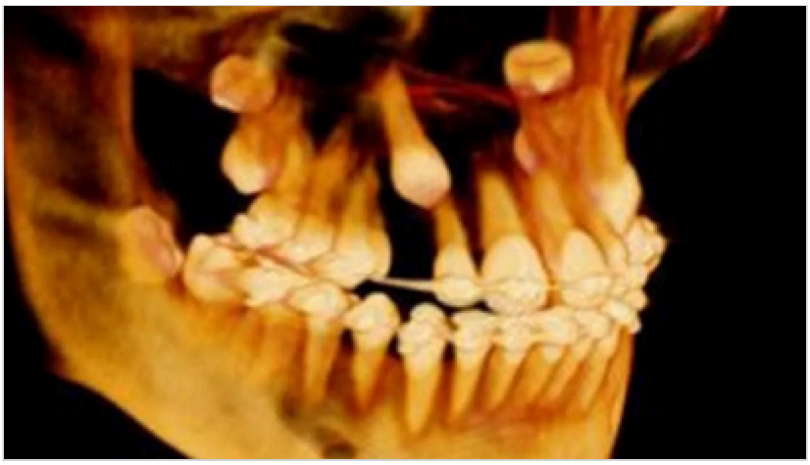


Endodontics


Using the i-CAT FLX V8, V10, or V17, clinicians can survey roots in three dimensions.
When more concentrated studies are necessary, high resolution scans – up to .125 voxels – lend more detail for the identification of lesions. Scans can also be collimated to cover the area of interest.

Within Tx STUDIO software, scans can be explored axially and buccolingual for a complete survey of fractures, accessory canals and endo-perio involvement.


Give Feedback




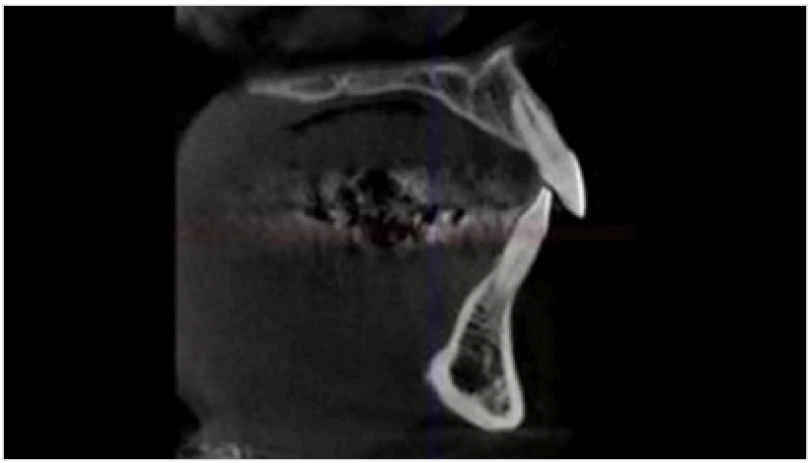

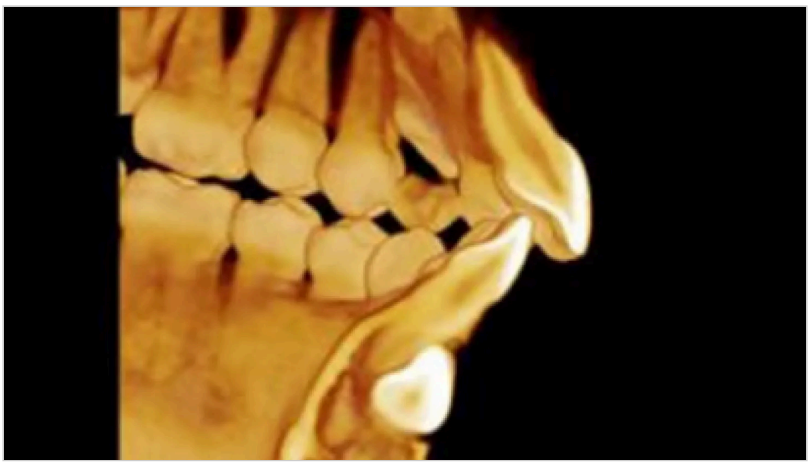

Give Feedback

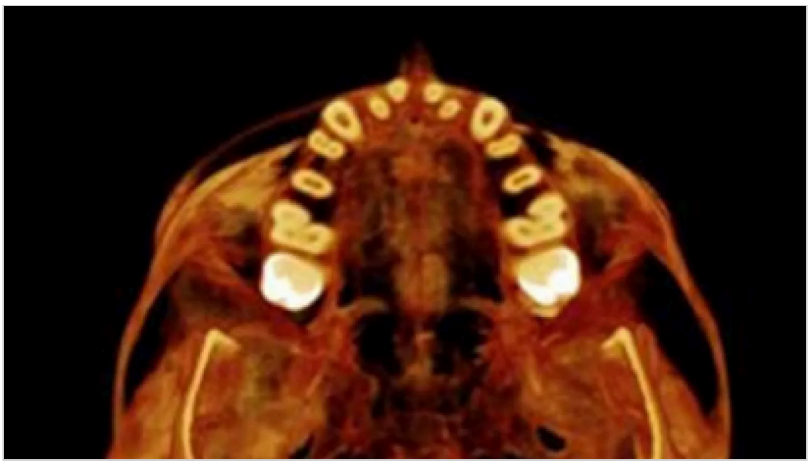

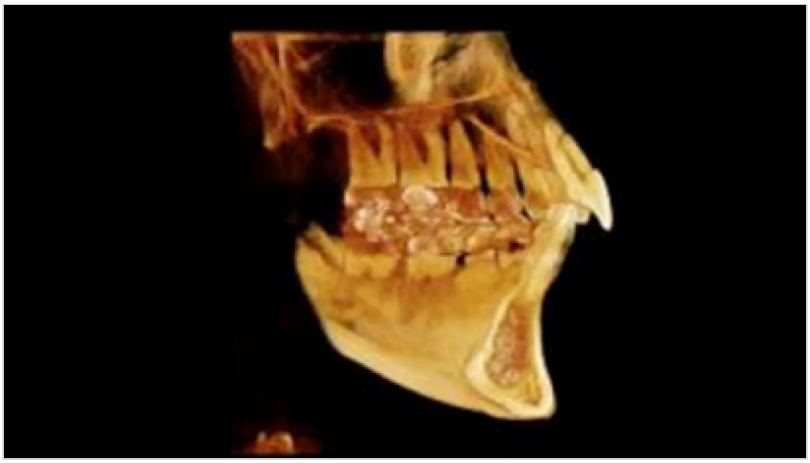


Give Feedback


Airway Analysis


Using the i-CAT FLX V10 or V17, clinicians can evaluate airway obstructions to help identify patients with an increased risk of a possible sleep disorder.

For patients with suspect airway or sinus tissues, you can use Tx STUDIO software to review the 3D data and to reveal restricted airways and determine appropriate treatments with precise anatomical views and measurements. Assess airway volume-at-a-glance using color-coded constriction values. Quickly trace airway on-screen to perform automatic calculations and measurements of paranasal sinuses to evaluate treatment options. i-CAT is FDA cleared to provide sufficient and appropriate data from 3D volumes of the maxillofacial anatomical structures, including the


airway.

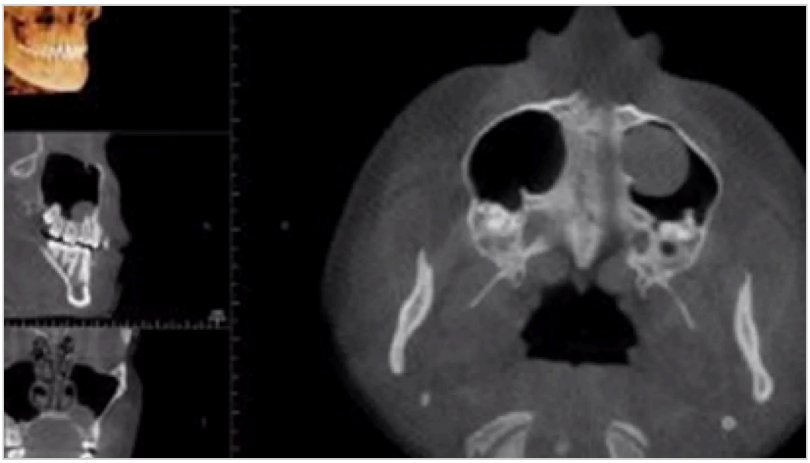


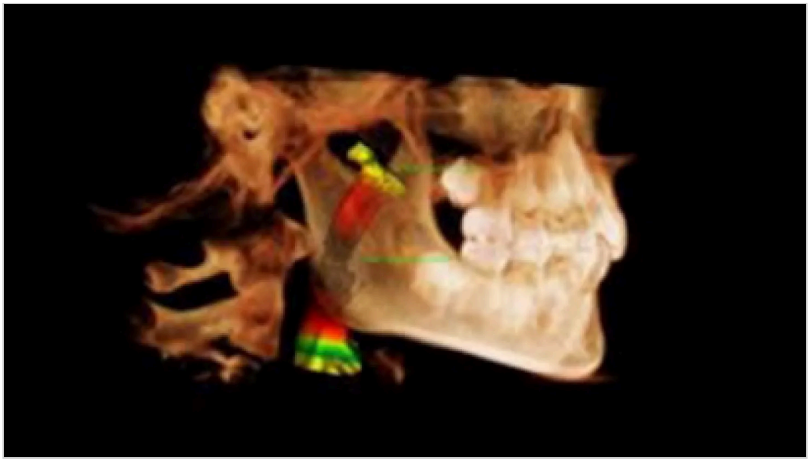


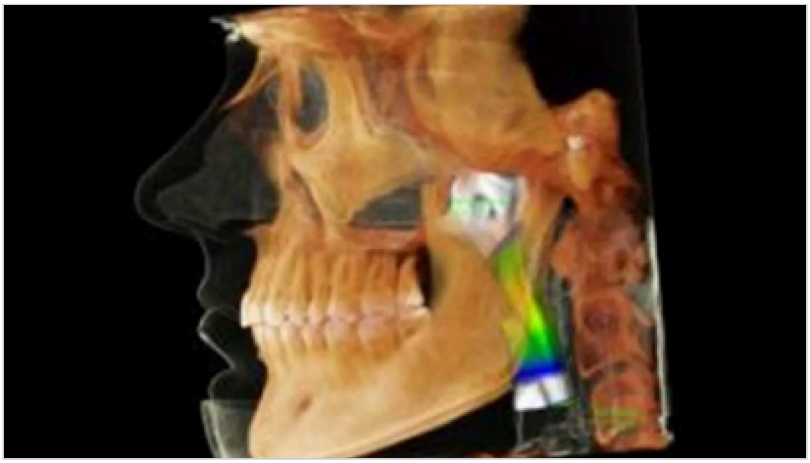


Give Feedback

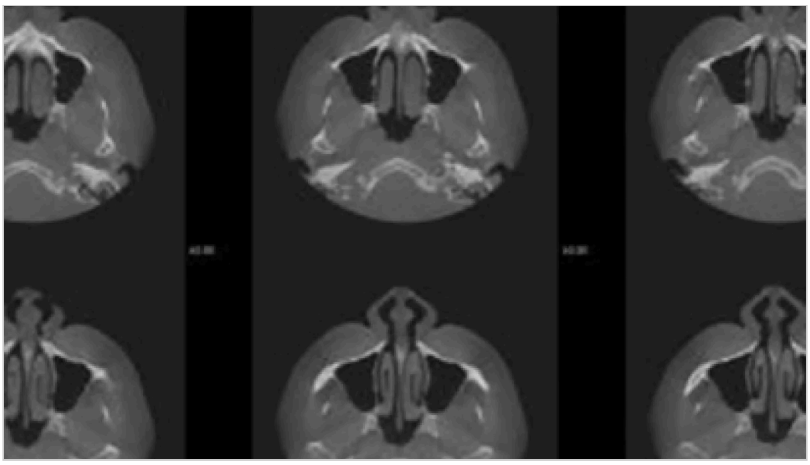

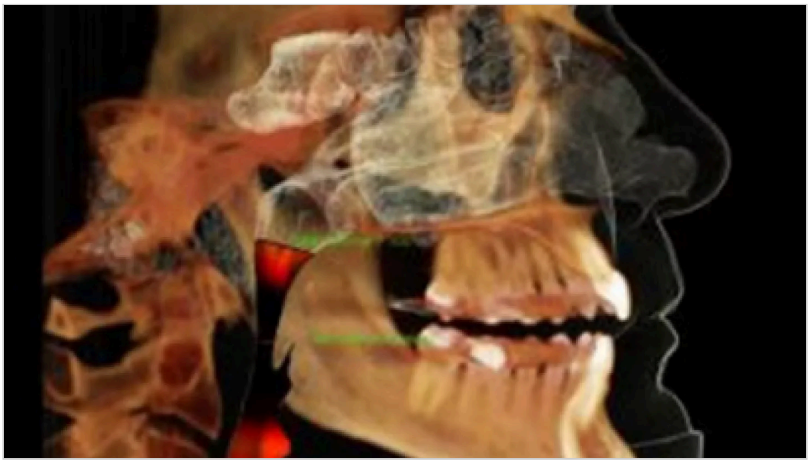

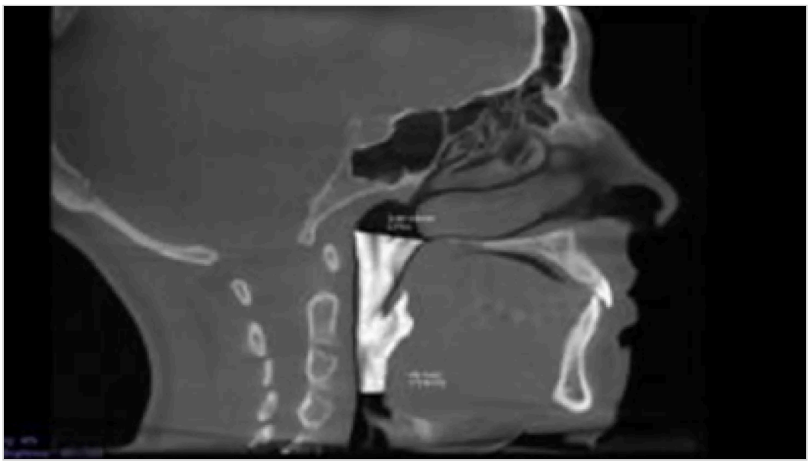

Give Feedback

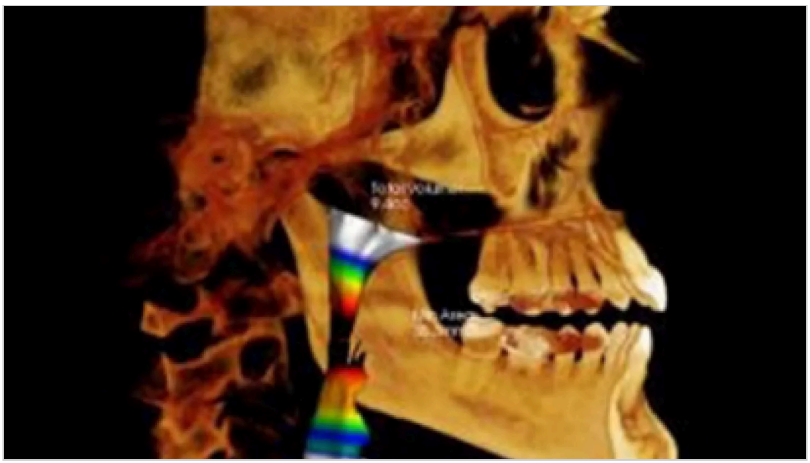



Periodontics


Using the i-CAT FLX V10 or V17, clinicians can analyze bone structure. From implant placement to surgical options for the management of bone loss, the i-CAT FLX provides periodontists with a wide range of services expected from the specialty. Capture 3D volumetric images for a more thorough analysis of bone structure as well as sinus and nerve location.

i-CAT imaging system provides options and a full complement of tools fo implant placement. Use scan data to help plan the course of treatment o bony defects prior to the actual osseous surgery appointment. With Tx STUDIO software's highly visual 3D presentation, share your diagnosis and plan with the patient for greater understanding leading to better post-treatment compliance.


Give Feedback

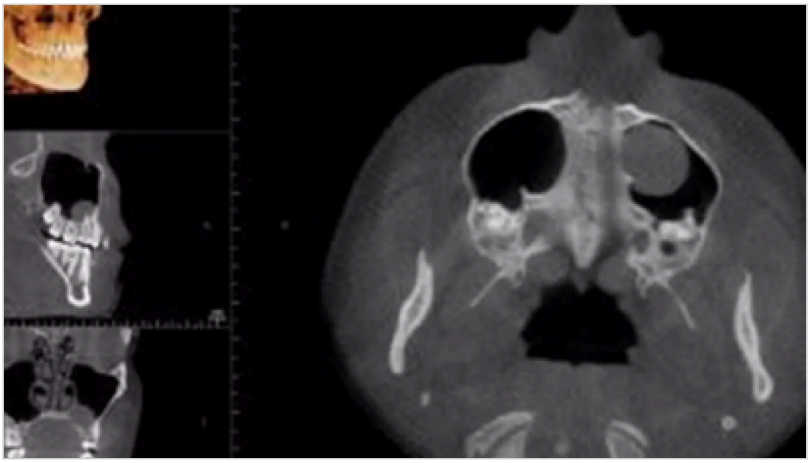


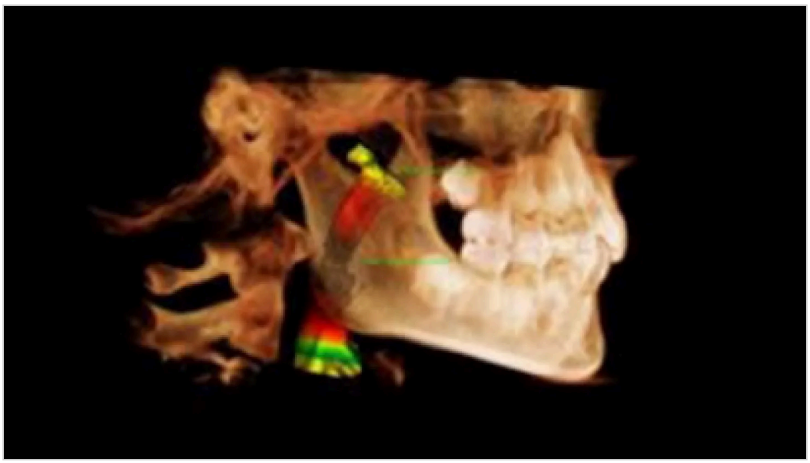


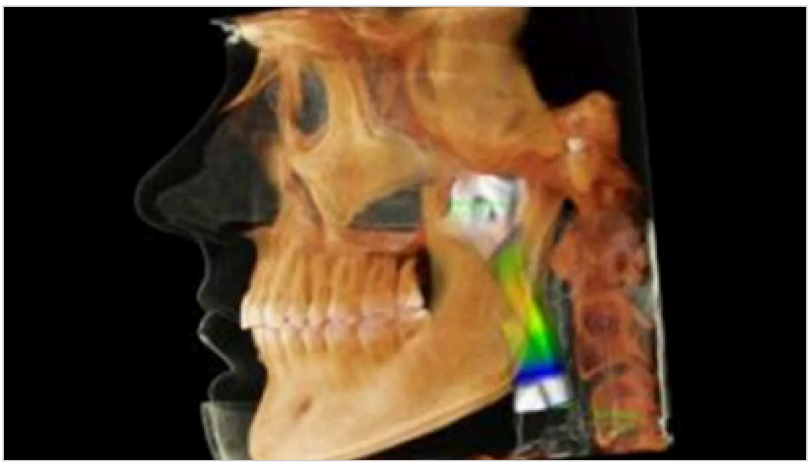


Give Feedback

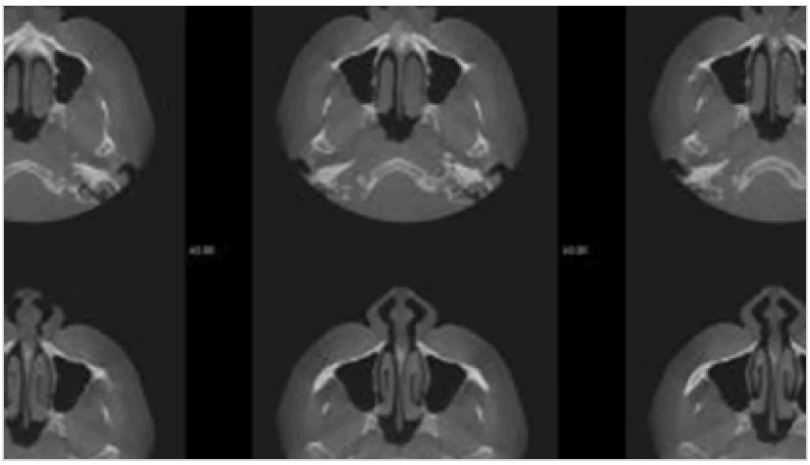

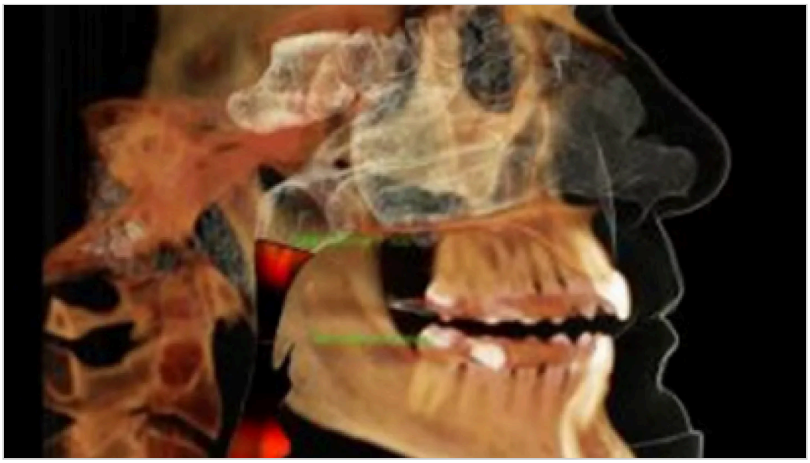

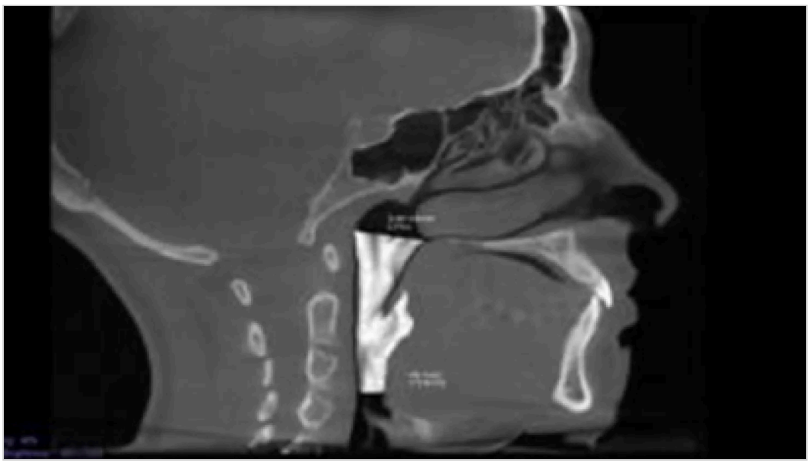

Give Feedback

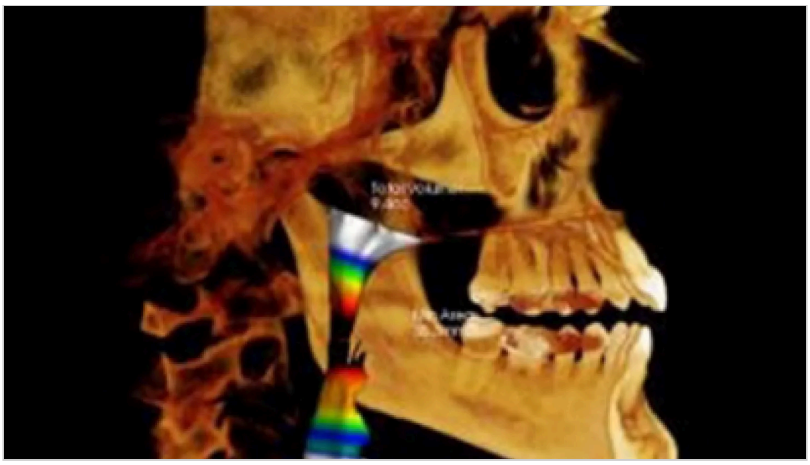



TMJ


Using the i-CAT FLX V10 or V17, clinicians can detect TMD and access fractures.Detect TMJ anomalies for the ability to design effective patient treatment. Using the TMJ visualization tools, zero in on the temporomandibular joints to identify wear, defects, and fractures. These tools also act as virtual study models to streamline and expedite treatment and allow you to design splints with the optional Medical Design Studio module.

Give Feedback




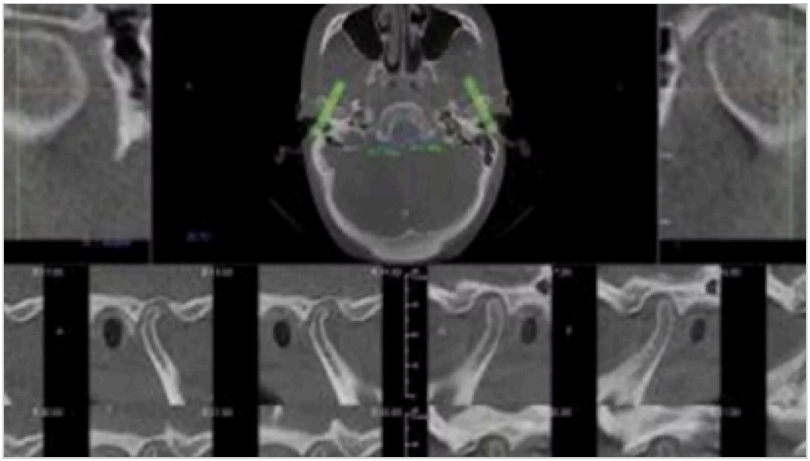

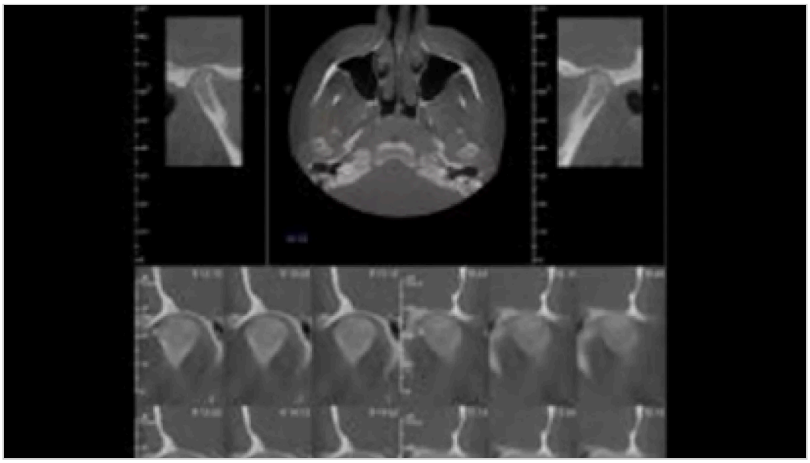

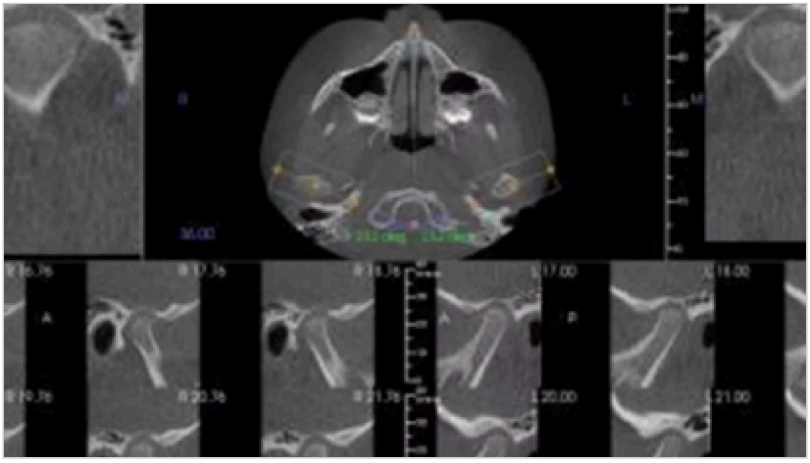

Give Feedback

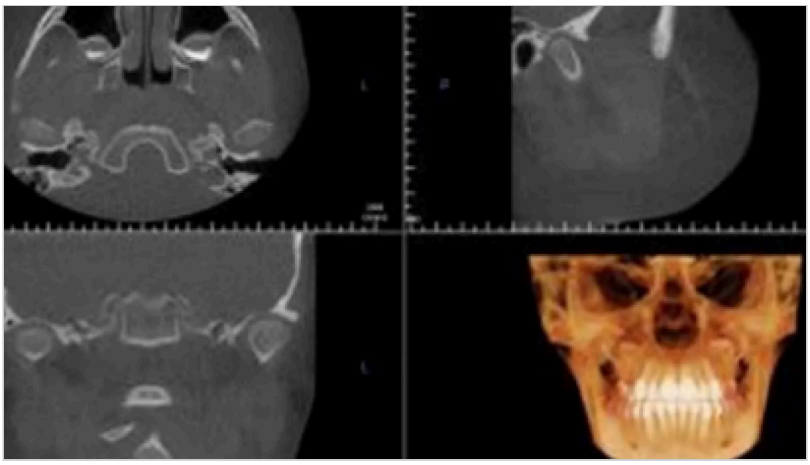





Give Feedback



Specifications

Give Feedback

i-CAT FLX Specifications

Sensor Type	Amorphous Silicon Flat Panel Sensor with CsI Scintillator
Grayscale resolution	16-bit
Voxel size	.4 mm, .3 mm, .25 mm, .2 mm, .125 mm
Collimation	Electronically controlled fully adjustable collimation
Scan time	4.8, 8.9, 14.7, 17.8 or 26.9 seconds*
Exposure Type	Pulsed
Field of view	V8: 8 cm x 5 cm, 8 cm x 8 cm V10: 8 cm x 5 cm, 8 cm x 8 cm, 16 cm x 4 cm, 16 cm x 6 cm, 16 cm x 8 cm, 16 cm x 10 cm V17: 8 cm x 5 cm, 8 cm x 8 cm, 16 cm x 4 cm, 16 cm x 6 cm, 16 cm x 8 cm, 16 cm x 10 cm,16 cm x 11 cm, 16 cm x 13 cm, 23 cm x 17 cm
Reconstruction shape	Cylinder
Typical reconstruction time	Less than 30 seconds
Viewing and treatment software	Included
DICOM compatible	Yes
Unit size	48” (w) x 69.5” (h) x 36.37” (d)
Patient position	Seated

*Utilizing the i-CAT FLX QuickScan+ exposure protocol. Use of lower dosage imaging may only be suitable for certain diagnostic tasks. Image quality is proportional to dose. i-CAT FLX offers a variety of exposure protocols allowing clinicians to adjust dosage to specific diagnostic needs.

X-ray images acquired using i-CAT can be analyzed within Tx Studio, powered by Anatomage, for treatment planning. NOTE: All clinical images shown are created from scans utilizing i-CAT technology.

Give Feedback

Videos



The i-CAT Story



Give Feedback

i-CAT Product Brochure

Downloads



i-CAT V-series Family Brochure

Document type: Product Brochure

Language: English (US)

Release date: Wednesday, August 2, 2017 - 12:00

File size: 15.98 MB

Download

Give Feedback

To view Documents, please visit the Download Center.

DX00889 RevA

PRODUCTS

COMPANY

SUPPORT

DEXIS™ solutions transform practices and patient smiles. Over 150,000 DEXIS dental imaging devices can be found in practices all over the world. Discover how DEXIS can deliver for you.

Subscribe to our newsletter

Subscribe

Get the latest industry news, product information and more.

Give Feedback

Join the conversation

Terms of Use Privacy Policy Cookie Policy Compliance & Transparency EULA Patents

United States

© Copyright Dental Imaging Technologies Corporation.
All Rights Reserved.