

EXHIBIT B-3

Invalidity of 11,644,693 in view of Digital Hearing Aid Systems

On information and belief, systems comprising a digital hearing aid, receiver, and remote transmitter (“Digital Hearing Aid systems” or “DHA systems”) were sold, on–sale, in public use, or otherwise publicly available prior to the earliest priority date of U.S. Patent No. 11,644,693 (“’693 patent”) and qualify as prior art to the ’693 patent at least under pre-AIA 35 U.S.C. §§ 102(a) and/or (b). DHA systems anticipate or, alone or with other references, render obvious, one or more of the Asserted Claims of the ’693 patent.

The chart below illustrates how exemplary DHA systems render the Asserted Claims of the ’693 patent invalid as anticipated or obvious over the DHA systems alone or in view of additional secondary references, including but not limited to the references listed below and prior art identified in Defendant’s Invalidity Contentions. One exemplary DHA system identified in this claim chart includes, for example, the Phonak Perseo digital hearing aid, MicroLink receiver, and Smartlink SX transmitter (“Phonak system”), which were sold, on–sale, in public use, or otherwise publicly available at least as early as March 2004. The Phonak system is evidenced by, for example, the Phonak Perseo Brochure, Phonak SmartLink SX User Guide, Phonak – SmartLink RX – FM Systems Brochure, and Phonak – MicroLink User Guide, which published at least as early as the earliest claimed priority date of the ’693 patent, as well as the additional information cited herein.

The citations and descriptions provided herein refer to the Phonak system as a representative system demonstrating the relevant prior art features and functionality of the DHA systems, and is not intended to limit the invalidity contentions to this specific implementation. On the contrary, as demonstrated herein, components of the DHA systems, including the receiver and transmitter, could be used interchangeably with different hearing aids or compatible systems. Additionally, the hearing aid device could have an integrated receiver within a unified assembly. Accordingly, the chart identifies other examples of exemplary components that could be substituted for or used interchangeably with the MicroLink receiver and/or Smartlink SX transmitter, including but not limited to the Oticon DigiFocus Compact digital hearing aid. Moreover, and as discussed in greater detail in the following chart, the prior art DHA systems may comprise one or more of the following components which could be used interchangeably and include substantially similar structures and functionality:

- Digital Hearing Aids: Phonak: Perseo, Perso with PersonalLogic, Claro 11, Claro 22, Claro 21dAZ, Claro dAZ, Claro dAZ FM, Supero 411, Supero 412, PicoForte3 SC-D, PicoForte3, SuperFront, Clasica, Aero 11, Aero 22, Aero 23 mZ, Aero 33, MicroZoom P2 AZ, Astro SC, Astro MC, Aero 211, Aero 211 AZ, Aero 311, Aero 311 AZ, Supero 412, Supero 413 AZ, Piconet 2 P2, Piconet P2 AZ, SonoForte P3 AZ, PowerZoom P4 AZ, NovoForte E3, NovoForte E4, Conversa, Unison; Oticon, Inc.: DigiFocus Compact, Syncro; Starkey Hearing Technologies: Genesis 4, Genesis DX with one or more of Precision Audio Imaging, Axent II, Sequel, Endeavour, Arista; Siemens Hearing Instruments, Inc.: Acuris with e2e wireless, Triano, Signia Select, Prisma 2, Music, Music Digital, Phoenix, Phoenix Pro, Infiniti 3; GN ReSound: ReSound Air; Sonic Innovations Inc.: Natura 3 with DIRECTIONALplus, Natura 2 SE, Altair, Tribute, Adesso; Widex: Senso Diva; Micro-Tech: Denali, Caballo, Dx3; Microsound A/S: Microsound Pilot with Eartunner; University Wisconsin + Nicolet Instrument Corp.: Project Phoenix; and any other prior art digital hearing aids devices, including, but not limited to devices offered by the following manufacturers: Oticon Philips, Phonak, Siemens, Sonic Innovations, and Widex.
- Receivers: Phonak: MicroLink (MLx S, ML 8, MLxi, ML3-ML7), MicroVox RX2
- Transmitters: Phonak: SmartLink, SmartLink SX, Campus S, ClearCall, SoundPilot FM, WatchPilot Classic, WatchPilot Elegance, MAXX, EduLink, CrosLink, HandyMic

On information and belief, other configurations, models, or embodiments may also satisfy one or more claim limitations, and the inclusion of the examples identified in this chart do not preclude the existence or applicability of alternative implementations. Defendant reserves the right to

supplement, amend, or otherwise update these contentions as additional information becomes available, including but not limited to information obtained through discovery.

In particular, Defendant may rely on information obtained through at least the following third-party subpoenas:

- Subpoena to Blackberry Corporation, dated May 14, 2026
- Subpoena to GN Audio USA Inc., dated May 14, 2026
- Subpoenas to International Business Machines Corp., dated May 14, 2026
- Subpoenas to WS Audiology USA, Inc., dated May 13, 2026
- Subpoenas to Sonova USA, Inc., dated May 13, 2026
- Subpoenas to Semtech Corp., dated May 13, 2026
- Subpoenas to Plantronics, Inc., dated May 13, 2026
- Subpoenas to Nuance Communications Inc., dated May 13, 2026
- Subpoenas to Oticon Inc., dated May 28, 2026

To the extent IngenioSpec contends the DHA systems do not disclose one or more limitations of the claims, it would have been obvious to combine the DHA systems with the knowledge of a person of ordinary skill in the art (“POSITA”) and/or with one or more of the following references to render the Asserted Claims the ’693 patent invalid, for example, for the same reasons and in the same ways as set forth in Exhibits B-1 and B-2:

- U.S. Patent No. 5,721,783 (“Anderson” or “Anderson ’783”) was filed on June 7, 1995, and issued on February 24, 1998, and is therefore prior art under 35 U.S.C. §§ 102 (a), (b), and (e).
- U.S. Patent No. 7,245,730 (“Enzmann” or “Enzmann ’730”) was filed on January 13, 2003, published on July 15, 2004, issued on July 17, 2007, and is thus prior art under at least 35 U.S.C. §§102(a) and (e). Enzmann ’730 issued from U.S. Application No. 10/341,293, which was published as U.S. Patent Publication No. 2004/0136555 (“Enzmann ’555”). Enzmann ’730 and Enzmann ’555 both qualify as prior art and have corresponding disclosures, and citations herein to Enzmann ’730 includes and incorporates by reference citations to the corresponding disclosures in Enzmann ’555.
- U.S. Patent No. 7,187,960 (“Abreu” or “Abreu ’960”) was filed on April 22, 2003, claims priority to a provisional application filed on April 22, 2002, and published on March 25, 2004, and is thus prior art under at least 35 U.S.C. §§102(a) and (e).
- U.S. Patent No. 7,620,433 (“Bodley” or “Bodley ’433”) was filed on January 15, 2003, published on July 15, 2004, and issued on November 17, 2009, and is thus prior art under at least 35 U.S.C. §§102(a) and (e).
- U.S. Patent No. 7,046,230 (“Zadesky” or “Zadesky ’230”) was filed on July 1, 2002, published on April 24, 2003, and issued on May 16, 2006, and is thus prior art under at least 35 U.S.C. §§102(a), (b), and (e).

- U.S. Patent No. 6,671,379 (“Nemirovski” or “Nemirovski ’379”) was filed on July 31, 2001, and issued on December 30, 2003, and is thus prior art under at least 35 U.S.C. §§102(a), (b), and (e).
- U.S. Patent Publication No. 2005/0058313 (“Victorian” or “Victorian ’313”) was filed on September 11, 2003, and published on March 17, 2005, and is thus prior art under at least 35 U.S.C. §§102(a) and (e).
- U.S. Patent No. 7,150,526 (“Jannard” or “Jannard ’526”) was filed on July 28, 2003 and is thus prior art under at least 35 U.S.C. §§102(a) and (e).
- WIPO Publication No. WO2001054458 (“Eaton” or “Eaton ’458”) was filed on January 18, 2001, and published on July 26, 2001, and is therefore prior art under 35 U.S.C. § 102(b).
- WIPO Publication No. 2004/093490A1 (“Hollemans” or “Hollemans ’490”) was filed on April 15, 2004, and claims priority to a PCT application filed on April 18, 2003, and is thus prior art under at least 35 U.S.C. §§102(a) and (e).
- U.S. Patent Publication No 2005/0136839A1 (“Seshadri ’839” or “Seshadri”) was filed on October 27, 2004, claims priority to a provisional application filed on May 28, 2003, and published on June 23, 2005, and is thus prior art under at least 35 U.S.C. §102(e).
- U.S. Patent Publication No 2005/0226446A1 (“Luo ’446” or “Luo”), filed on April 8, 2004 and published on October 13, 2005, and is thus prior art under at least 35 U.S.C. §102(e).
- U.S. Patent No. 7010332B1 (“Irvin ’332” or “Irvin”), filed on February 21, 2000, was patented on March 7, 2006, and is thus prior art under at least 35 U.S.C. §102(e).
- U.S. Pat. No. 6,373,942 (“Braund”) was filed on April 7, 2000 and issued on April 16, 2002, and is thus prior art under at least 35 U.S.C. §§ 102 (a), (b), and (e).
- U.S. Pat. No. 7,979,035 (“Griffin”) was filed on Dec. 4, 2007 with a priority claim to Nov. 7, 2000, and is thus prior art under at least 35 U.S.C. § 102 (e).
- U.S. Pub. No. 2003/0068057 (“Miller”) was filed on Oct. 6, 2001 and published on Apr. 10, 2003, and is thus prior art under at least 35 U.S.C. §§ 102 (a), (b), and (e).
- U.S. Pat. No. 4,955,729 (“Marx”) was filed on Mar. 29, 1988 and issued on Sep. 11, 1990, and is thus prior art under at least 35 U.S.C. §§ 102 (a), (b), and (e).
- U.S. Pat. No. 7,542,012 (“Kato”) was filed on June 16, 2005 with a priority claim to Dec. 22, 2003 and issued on Jun. 2, 2009, and is thus prior art under at least 35 U.S.C. § 102(e).

For any claim limitation, Defendant may rely on excerpts cited for any other limitation and/or additional excerpts not set forth fully herein to the extent necessary to provide further explanation for a reference’s disclosure of a limitation. Where an excerpt refers to or discusses a figure or figure items, that figure and any additional descriptions of that figure should be understood to be incorporated by reference as if set forth fully herein.

With respect to the proposed modifications to the DHA systems, such modification would have been obvious to try, an obvious combination of prior art elements according to known methods to yield predictable results, a simple substitution of one known element for another to obtain predictable results, a use of known techniques to improve a similar devices or method in the same way, an application of a known technique to a known device or method ready for improvement to yield predictable results, a variation of a known work in one field of endeavor for use in either the same field or

a different one based on design incentives or other market forces with variations that are predictable to a POSITA, and/or obvious in view of teachings, suggestions, and motivations in the prior art that would have led a POSITA to modify or combine the prior art references.

Defendant expressly reserves the right to rely on additional prior art references to demonstrate that certain features were well known in the art and would have been obvious to a POSITA as of the effective filing date of the '693 patent, e.g., to the extent the references cited in the accompanying chart are argued or found not to expressly teach or suggest such features. For example, use of directional microphones in headworn devices, e.g., for voice recognition, receipt of audio commands, sensing of environmental sounds, and noise cancellation were well known features, and it would have been obvious to a POSITA to combine them with those of the DHA systems for the reasons set forth below. For example, Braund discloses well-known aspects of using a microphone for voice recognition (see, e.g., Braund at 3:47-58, 6:23-27, 8:55-61) and rechargeable battery features (see, e.g., Braund at 9:6-7), as recited in at least claims 17, 25, 26, 52, 83, and 72. Griffin discloses well-known aspects of directional microphones for noise cancellation (see, e.g., Griffin at 4:21-39) and rechargeable battery features (see, e.g., 4:21-5:8, 9:26-33), as recited in at least claims 6-12, 17, 50, 51, 53, 61, and 64-72, 25, 26, 52, and 72. Miller discloses well-known aspects of directional microphones for noise cancellation (see, e.g., Miller at paras. [0044]-[0045], [0062]), as recited in at least claims 6-12, 17, 50, 51, 53, 61, and 64-72. Other well-known features include use of sensors in, on, or with headworn devices, e.g., to determine if/when headworn devices were being worn. For example, Marx discloses well-known aspects of using sensor data to determine if a device is being worn and modifying operation of the device based on that determination (see, e.g., 3:34-4:12), as recited in at least claims 23, 54, 55, 56, 57, and 73. For further example, Kato discloses well-known aspects of including an accelerometer in a head worn device (see, e.g., 35:3-9), as recited in at least claim 30.

To the extent that these Invalidity Contentions rely on or otherwise embody particular interpretations of terms or phrases in the asserted claims, Defendant is not proposing any such interpretations as proper interpretations or constructions of those terms or phrases. Various positions put forth in this document are predicated on IngenioSpec’s apparent interpretation of its claims as evidenced by its Infringement Contentions. Those positions are not intended to and do not necessarily reflect Defendant’s interpretation of the true and proper scope of the asserted claims. Defendant further incorporates by reference all reservations of rights set forth in the Invalidity Contentions, as if fully set forth herein.

Defendant’s arguments throughout these invalidity contentions are based on IngenioSpec’s apparent interpretation of the claims in its infringement contentions. Defendants do not agree that such interpretation is proper and reserve the right to argue that the mapped actions do not infringe any claims of the '693 patent. In the alternative, should the Court or factfinder construe and/or apply the claims in a manner consistent with IngenioSpec’s interpretations in its infringement contentions for the '693 patent, the claims are invalid for at least the reasons set forth herein.

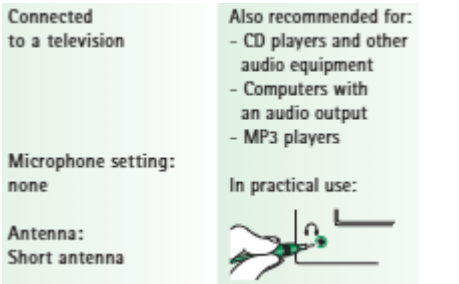
Claims	Prior Art
1[pre]. A personal audio system operable to at least generate customized audio output, comprising:	DHA systems, alone or in view of Anderson or Enzmann, expressly or inherently discloses “[a] personal audio system operable to at least generate customized audio output.”

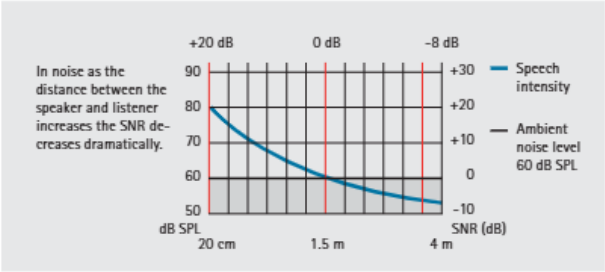
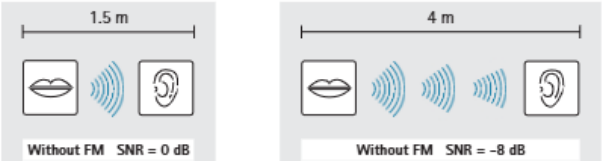
Claims	Prior Art
	<p data-bbox="583 175 1486 207">Phonak - Hearing Aid Museum re SmartLink SX FM System at Page 1</p> <p data-bbox="638 232 1050 410">The SmartLink SX was a three-in-one device: an FM transmitter, a mobile phone bluetooth link and hearing aid remote control. The SmartLink SX worked with all Phonak's hearing aids beginning with the Perseo, Claro and Supero.</p> <p data-bbox="583 440 1268 472">Phonak - SmartLink SX - FM Systems Manual at p. 3</p> <p data-bbox="638 492 1092 792">SmartLink, based on Phonak breakthrough digital FM synthesizer technology, now enables people with a hearing impairment and those with normal hearing to use the same communication tools with equal effectiveness. SmartLink is a communication gateway to all situations where understanding perfectly, all of the time, is key. The unique innovation, SmartLink, meets the needs of people with hearing loss, giving them access to effective communication in difficult situations and to the many advantages off modern communication tools, via one device.</p>

Claims	Prior Art
	<p data-bbox="583 175 1272 211">Phonak - SmartLink SX - FM Systems Manual at p. 4</p>  <p>The diagram illustrates a mobile phone with three input/output ports on its side. Arrows indicate connections: 'Speech' (a person's mouth) connects to the top port, 'Bluetooth link' (a small device) connects to the middle port, and 'Audio input' (a microphone) connects to the bottom port. The phone is shown in a vertical orientation with a silver finish.</p>

Claims	Prior Art
	<p data-bbox="583 175 1289 214">Phonak - SmartLink SX - FM Systems Manual at p. 10</p>  <p data-bbox="630 795 651 828">0</p> <p>The diagram illustrates the Phonak SmartLink SX FM system. It features a base unit (right) and a handset (left). The base unit has a screen displaying '22' and 'FM', and buttons labeled '+', 'P', '-', and 'F'. The handset has a screen displaying a menu with options: 'Exit smartlink mode', 'Move', 'Call class', 'Help', 'Exit', 'Select', and 'Cancel'. Arrows indicate the following connections: 'Bluetooth' between the handset and base unit; 'FM Signal' from the base unit to a user's ear; and 'Speech' from the user's mouth to the base unit. Labels 'Press to accept call' and 'Press to end call' point to the '+' and 'P' buttons on the base unit, respectively.</p>

Claims	Prior Art
	<p data-bbox="583 175 1289 207">Phonak - SmartLink SX - FM Systems Manual at p. 14</p> <div data-bbox="627 224 1285 823">  <p data-bbox="667 446 1186 527">SmartLink™ SX</p> <p data-bbox="667 552 1029 584">The gateway to communication</p> <ul data-bbox="667 649 1081 812" style="list-style-type: none"> • The only FM system with digital synthesizer • Cellular phone access via Bluetooth • Integrated remote control for hearing instruments • Ergonomic design with visual controls • State-of-the-art energy management </div> <p data-bbox="583 841 1136 873">Phonak - SmartLink SX User Guide at p. 6</p> <p data-bbox="636 893 1062 966">This unique combination of an FM transmitter, mobile phone link and hearing instrument remote control is the gateway to a completely new hearing experience.</p>

Claims	Prior Art
	<p data-bbox="583 175 1150 207">Phonak - SmartLink SX User Guide at p. 32</p>  <p data-bbox="583 672 1150 704">Phonak - SmartLink SX User Guide at p. 47</p>  <p data-bbox="583 1024 1150 1057">Phonak - SmartLink SX User Guide at p. 58</p> <p data-bbox="638 1078 1052 1149">The SmartLink is used by hearing-impaired people. To operate the SmartLink, BTE hearing instruments with MicroLink receivers are required.</p>

Claims	Prior Art
	<p>Phonak – Perseo Brochure at Page 5</p> <p>The Science of PersonalLogic™</p> <p>Phonak has applied knowledge from new research and past experience to add a new dimension of individualized logic to automatic hearing systems – PersonalLogic. Combining new and proven features, PersonalLogic optimizes the functioning, operation and fitting of hearing systems for each individual user.</p> <p>Thanks to PersonalLogic, Perseo takes customization to a level beyond "custom-fitted" related to individual hearing loss. The new digital system from Phonak is the first hearing system with customized automatic operation, as opposed to operation based on factory set, unchangeable rules. Now, Perseo gives hearing professionals control over automatic operation, allowing them to provide clients with a hearing system "custom-tailored" to their personal hearing requirements and tastes.</p> <p>Phonak – Perseo Brochure at Page 15</p> <p>The perfectly integrated communication system</p> <p>FM MicroLink</p> <p>FM technology has allowed thousands of children to participate in standard education settings and to be actively involved in social activities. Phonak pioneered this development and the integration of FM into hearing instruments to produce communication systems that are effective in even the most demanding situations. Since the innovation of miniaturized Phonak MicroLink systems, FM has become an important solution also for adults with a high communication need in work and leisure environments.</p> <p>All Perseo BTEs have been designed to incorporate MicroLink with its fully integrated appearance and operation. This integrated system is unique to Phonak.</p> <p>Speech-to-Noise Ratio (SNR) as a function of distance between speaker and listener</p>  <p>In noise as the distance between the speaker and listener increases the SNR decreases dramatically.</p> <p>Speech-to-Noise Ratio without FM is distance dependent</p>  <p>Without FM SNR = 0 dB</p> <p>Without FM SNR = -8 dB</p>

Claims

Prior Art

Phonak – Perseo Brochure at Page 26

Perseo 11

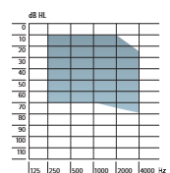


||| Digital microphone

OPTIONS
||| Tactronic switch
||| Control beeps

Ear Simulator data
45 dB max. gain
119 dB SPL max. output

2 cc data
35 dB max. gain
108 dB SPL max. output



Perseo 12

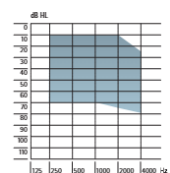


||| Digital microphone
||| Volume control

OPTIONS
||| Tactronic switch
||| Control beeps

Ear Simulator data
50 dB max. gain
119 dB SPL max. output

2 cc data
40 dB max. gain
108 dB SPL max. output



Perseo 22

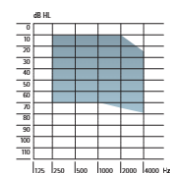


||| Telecoil

OPTIONS
||| Tactronic switch
||| Control beeps
||| Remote control WatchPlot2 and SoundPlot

Ear Simulator data
50 dB max. gain
119 dB SPL max. output

2 cc data
40 dB max. gain
108 dB SPL max. output



Claims

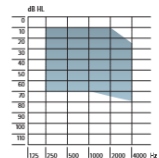
Prior Art

Phonak – Perseo Brochure at Page 27

Perseo 23 dAZ



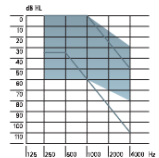
- Adaptive digital AudioZoom
- Telecoil
- OPTIONS
 - Tactronic switch
 - Control beeps
 - Remote control WatchPlot2 and SoundPilot
- Ear Simulator data
 - 58 dB max. gain
 - 118 dB SPL max. output
- 2 cc data
 - 45 dB max. gain
 - 108 dB SPL max. output



Perseo 111 dAZ
Perseo 111 dAZ FM



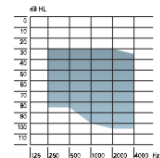
- FM solutions
- Adaptive digital AudioZoom
- Telecoil
- Tactronic switch
- OPTIONS
 - Control beeps
 - Remote control WatchPlot2 and SoundPilot
- Ear Simulator data
 - 54 dB max. gain
 - 127 dB SPL max. output
- 2 cc data
 - 47 dB max. gain
 - 121 dB SPL max. output



Perseo 211 dAZ
Perseo 211 dAZ FM



- FM solutions
- Adaptive digital AudioZoom
- Telecoil
- Tactronic switch
- OPTIONS
 - Control beeps
 - Remote control WatchPlot2 and SoundPilot
- Ear Simulator data
 - 62 dB max. gain
 - 130 dB SPL max. output
- 2 cc data
 - 55 dB max. gain
 - 125 dB SPL max. output



Technical Data

Claims

Prior Art

Phonak – Perseo Brochure at Page 28

Perseo 311 dAZ Forte
Perseo 311 dAZ Forte FM



Remote controls



FM Transmitters



- FM solutions
- Adaptive digital AudioZoom
- Telecoil
- Tactronic Switch

OPTIONS

- Control beeps
- Remote control WatchPilot2 and SoundPilot

Ear Simulator data

72 dB max. gain
135 dB SPL max. output

2 cc data*

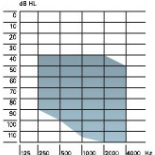
85 dB max. gain
130 dB SPL max. output

WatchPilot™2



HandyMic™





TelCom™



SoundPilot™



Campus S™



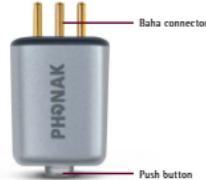
Phonak – Microlink User Guide at p. 24

MicroLink MLx can be used with Behind-The-Ear (BTE) hearing instruments from several different manufacturers.

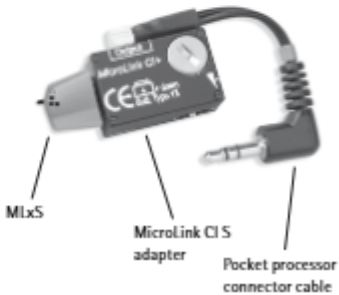
Phonak – Microlink User Guide at p. 52

The HandyMic and MicroVox TX2 are used by hearing impaired people. To operate the HandyMic and MicroVox TX2, hearing instruments with FM receivers (e.g., MicroLink) are required. The user benefit of FM is an improved signal-to-noise ratio which leads to an optimized speech understanding even in very difficult listening situations.

Claims	Prior Art
	<p>Phonak – Annual Report (2003–2004) at p. 12</p> <p>In addition to the very successful launch of Phonak’s new premium product line Perseo in March 2003, the following is a summary of the key product launches during the past fiscal year:</p> <ul style="list-style-type: none"> • MAXX and Amio, Phonak’s new entry level digital products, launched in June 2003; • Conversa, Unitron Hearing’s mid-level digital product line launched in June 2003; • Supero +, an enhanced version of Phonak’s digital high power line, launched in October 2003; • MAXX directional, which adds directional signal processing to the entry-level product family, launched in October 2003; and • Liaison, Unitron Hearing’s new high-end product line, launched in October 2003. <p>Phonak – Annual Report (2003–2004) at p. 12</p> <p>Within the FM business segment, the Group continued to leverage existing key technology into new applications. With SmartLink, Phonak has achieved another industry first by allowing hearing instrument users to access other communications systems. SmartLink combines a remote control, a transmitter, a directional microphone system and a Bluetooth link to mobile phone all in one single device. Phonak also presented EduLink, an innovative FM device to assist the learning development of children with Auditory Processing Disorder (APD). EduLink was perceived very positively within the medical community and recent classroom tests show favourable results.</p> <p>Phonak – Online Interview re Perseo (April 7, 2003) at p. .5</p> <p><u>AO/Beck</u>: Phonak has been a leader in FM technology-is this available for the Perseo product line?</p> <p><u>Dave Fabry/Phonak</u>: Yes, FM can be an integrated solution for the Perseo BTEs and is compatible with all Phonak transmitters.</p>

Claims	Prior Art
	<p>Phonak – Online Interview re Perseo (April 7, 2003) at Page 5</p> <p>In a nutshell, Perseo is designed for new and experienced users who are looking for a hearing instrument with advanced adaptive features, enhanced sound quality, easy personalization of fitting and choice of user operation (namely, automatic, manual program control with Tactronic switch, or via remote control). The choice of one product over another is dependent on the patient's hearing loss, listening needs, and prioritization of size, cost, and available performance features. Phonak provides an entire family of digital options from Perseo, Claro, Aero, and our new digital value line that will be introduced this summer. In addition, Supero is the only digital hearing instrument designed specifically for super-power applications.</p> <p>Phonak – Fitting Guide for MLxi and MLxi Baha at Page 2</p> <div data-bbox="636 609 1291 1209"> <div data-bbox="636 609 976 1209"> <p>Welcome</p> <p>MLxi is a universal Dynamic FM receiver. Compatible with virtually all BHE hearing instrument brands and types, it offers unparalleled ease of operation and is simply a case of plug and listen. Featuring totally automatic DPAI management, it also reduces the need for electro acoustic verification.</p> <p>Designed to offer the ultimate in performance, MLxi Baha has the same functionality as MLxi, but incorporates a small connector for Cochlear's Intenso, Divino and Compact Baha instruments.</p> <p>The state-of-the-art, highly miniaturized technology inside MLxi/MLxi Baha enables users to benefit from features never seen before:</p> <ul style="list-style-type: none"> • Dynamic FM • AutoConnect • Monitoring <p>Please read this fitting guide carefully to make sure the MLxi/MLxi Baha is well fitted to your customers hearing system so they can benefit from all its features.</p> <p>The following symbol appears in this user guide:</p> <p> Important information for the fitter.</p> <p>For detailed information on using the FM transmitter, which is part of the overall FM solution, please refer to the appropriate FM transmitter user guide (also see www.phonak.com).</p> <p>MLxi/MLxi Baha are quality products developed by Phonak, the world leader in FM technology, innovation and reliability.</p> </div> <div data-bbox="976 609 1291 1209"> <p>2. Instrument description</p> <p>2.1 MLxi</p>  <p>2.2 MLxi Baha</p>  </div> </div>

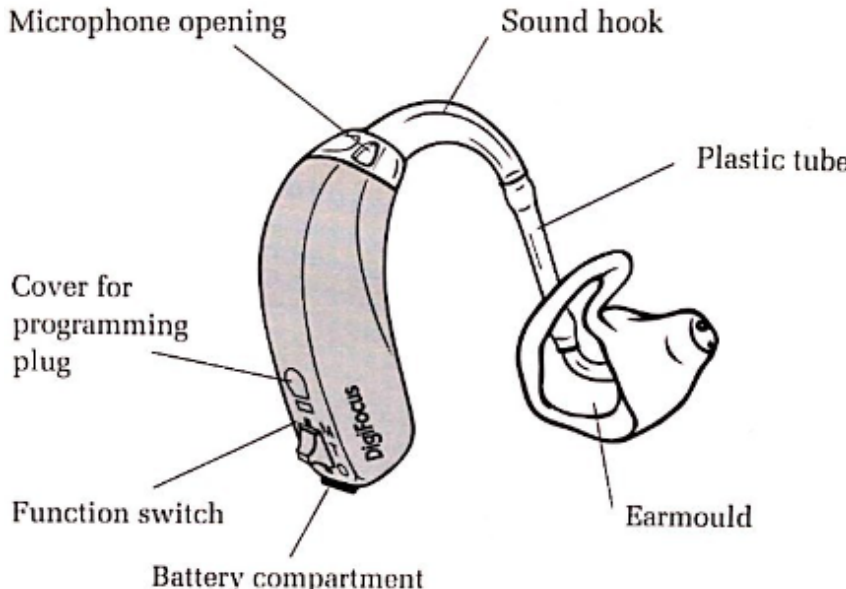
Claims	Prior Art
	<p>Phonak – Fitting Guide for MLxi and MLxi Baha at Page 9</p> <hr/> <p>Dynamic FM transmitters</p> <div data-bbox="636 267 1570 602"> <div>  <p>inspiro is a Dynamic FM transmitter for teachers and suitable for use in both regular and special schools.</p> </div> <div>  <p>DynaMic is the passaround microphone for the MultiTalker Network with inspiro.</p> </div> <div>  <p>SmartLink+ is the FM transmitter for those who want nothing less than the very best.</p> </div> <div>  <p>ZoomLink+ is the FM transmitter for those who know what they want.</p> </div> <div>  <p>EasyLink+ is the FM transmitter that is easiest to operate.</p> </div> </div> <p>Phonak – Claro 211 Product Information at Page 2</p> <hr/> <p>Principle of FM operation</p> <p>The transmitter unit is worn by the person speaking, placed in front of the sound source or even connected directly to TV or HiFi. The transmitter unit collects the desired signal and transmits it via radio waves directly to the FM module integrated in the hearing instrument. The incoming FM signal received by the hearing instrument is processed either in accordance with the audiological setting of the instrument, providing consistent sound quality or according to an individual FM program.</p>

Claims	Prior Art
	<p data-bbox="583 175 1226 207">Phonak SmartLink SX Product Information at p. 1</p> <div data-bbox="630 224 1129 818"> <p data-bbox="642 230 856 256">Product Information</p> <p data-bbox="642 263 982 321">SmartLink™ SX</p>  <p data-bbox="642 659 1121 815">SmartLink SX is a compact communication device for people with hearing loss who lead busy, active lives. The SmartLink SX incorporates an FM transmitter with digital signal processing, a remote control for Phonak hearing instruments, and Bluetooth for communication with mobile phones. With its digital frequency synthesizer, the SmartLink allows a highly flexible choice of operating frequency. For maximum customer benefit, we recommend using a Phonak multi-channel FM receiver (e.g. the MLxS or MLBS) that can be automatically synchronized from the SmartLink.</p> </div> <p data-bbox="583 841 1306 873">Phonak MicroLink MLxS Product Information at Page 1</p> <p data-bbox="638 896 1495 977">The MicroLink MLxS is the first miniaturized FM receiver with a digital frequency synthesizer for full multi-frequency flexibility. The receiver features fully automatic, non-contact frequency synchronization with either a WallPilot or the Campus S.</p> <p data-bbox="638 997 1453 1049">The MicroLink MLxS is compatible with all Phonak FM transmitters and the WallPilot, and is recommended for use with this equipment.</p> <div data-bbox="638 1107 1478 1451"> <div data-bbox="638 1107 945 1419"> <p data-bbox="638 1107 840 1127">MLxS with BTE instrument</p>  <p data-bbox="827 1237 945 1279">Audio shoe (e.g. AS4-MLxS)</p> <p data-bbox="869 1312 911 1331">MLxS</p> </div> <div data-bbox="1142 1107 1478 1451"> <p data-bbox="1142 1107 1453 1127">MLxS with adapter for CI pocket processor</p>  <p data-bbox="1142 1360 1184 1380">MLxS</p> <p data-bbox="1268 1377 1373 1419">MicroLink CI S adapter</p> <p data-bbox="1360 1416 1478 1451">Pocket processor connector cable</p> </div> </div>

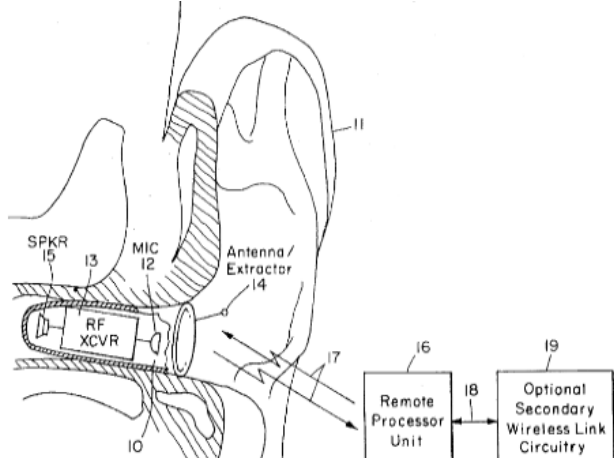
Claims	Prior Art
	<p>Phonak MicroLink MLxS Product Information at Page. 1</p> <p>• Compatible with most behind-the-ear hearing aids on the market</p> <p>Phonak – Collections 2002 Product List at Page 29</p> <div data-bbox="632 332 1071 930"> <p>Phonak's FM solutions: For maximum speech intelligibility</p> <div> <p>Whenever even latest-generation hearing instruments approach their physical limits – noise, reverberation, or distance – count on Phonak's FM solutions to deliver maximum speech intelligibility.</p> <p>The frequent difficulties encountered by hearing instrument users when watching TV and making phone calls can be overcome with Phonak's FM technology.</p> </div> <div> <p>Phonak FM solutions improve the SNR by 12 to 20 dB.</p> <p>The FM systems by Phonak can be used with the majority of commercially available BTE instruments and cochlear implants.</p> </div> <p>Recommended applications:</p> <div> <p> Dialog: For example at home, on a busy street, at the shopping mall or in a restaurant.</p> <p> Group discussions: For example at a family get-together, in a very noisy restaurant, at business meetings or conferences.</p> <p> Making telephone calls: The Phonak FM system provides access to the many capabilities of modern telecommunications devices.</p> <p> Watching TV and listening to music: Phonak FM provides access to the entire spectrum of today's information and entertainment technologies.</p> </div> <div> <p> Schools and continuing education: For example in classes at schools and universities, on outings and excursions, in lectures or on museum tours.</p> <p> Transportation: For example in cars, trains, trams and buses.</p> <p> Sports of all kinds: For example when cycling, playing golf, skiing, jogging or hiking.</p> <p> Events: For example at the theater, in concerts, at the movies, in church services or at community meetings.</p> </div> </div>

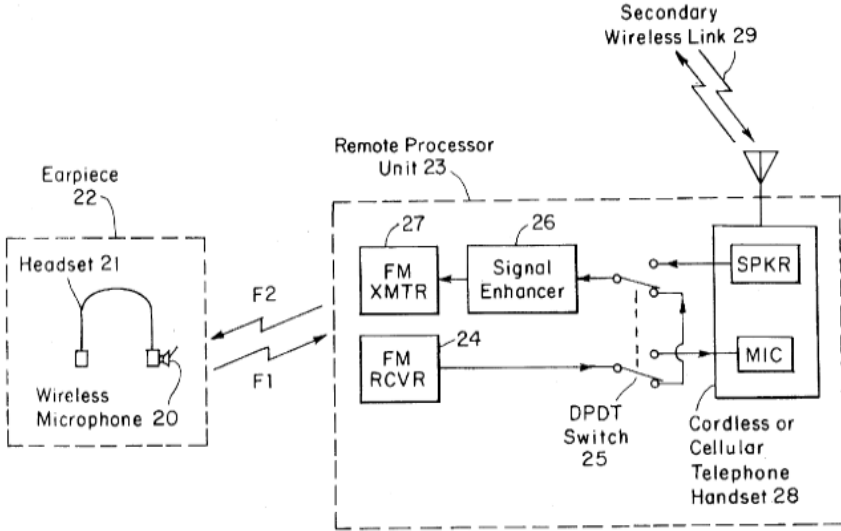
Claims	Prior Art
	<p>Phonak – Collections 2002 Product List at Page 30</p> <div><div><p>MLx compatible BTE makes:</p><ul style="list-style-type: none">ArgosyAudifonBeltoneDanavoxGN ReSoundHansatonIntertonOticonPhilipsPhonakRextonSiemensSonic InnovationsUnitronWidex</div><div><p>MLx, the most flexible FM receiver ever made, can be used not only with Phonak BTEs but also with most third-party hearing instruments available on the market today. This is possible thanks to the standard Euro connector which simply plugs into the audio shoe of the hearing instrument.</p><p>Once again, Phonak demonstrates its commitment to innovative developments that serve hearing-impaired adults and children.</p><p>Detailed, up-to-date compatibility lists of Phonak's MLx MicroLink receiver are available online at: www.phonak.com/mlx</p></div><div></div></div>

Claims	Prior Art
	<p data-bbox="581 175 1610 212">Lewis – Speech perception in noise: directional microphones (2004) at Page 431</p> <p data-bbox="636 228 1016 818">All subjects were fit with digital Phonak Claro 311 dAZ BTE hearing aids bilaterally. All earmolds had select-a-vent (SAV) venting and #13 or 3 mm horn tubing. In addition to the hearing aids, subjects were fit with Phonak Microlink ML8 FM receivers bilaterally. These FM receivers attach to the bottom of a BTE hearing aid and may be utilized in either the “FM only” mode, which attenuates the hearing aid microphone by 20 dB, or in the FM plus hearing aid mode, which allows for FM input and input of environmental sounds via the hearing aid microphones simultaneously at the same output level. The Phonak TX3 HandyMic FM transmitter served as the FM transmitter. This transmitter has three microphone options: (1) Wide Angle, which picks up sounds arriving from all directions around the transmitter microphone equally; (2) Zoom, which provides reduced amplification to signals arriving from the rear (cardioid); and (3) SuperZoom, which provides reduced amplification to signals originating from the rear and the sides (hypercardioid). The</p>

Claims	Prior Art
	<p data-bbox="583 175 1449 207">Oticon – DigiFocus Compact User Instructional Brochure at Page 4</p>  <p data-bbox="583 841 1079 873">Oticon – DigiFocus Description at p. 1</p> <p data-bbox="634 896 1314 967">DigiFocus is available both as behind-the-ear models (BTE) in colours which match the hair tone, and as in-the-ear (ITE) models which fit discreetly into the ear.</p>  <p data-bbox="714 1289 1239 1315"><i>DigiFocus Compact, Canal, Half Shell and Low Profile.</i></p>

Claims	Prior Art
	<p>Oticon – DigiFocus White Paper (1996) at Page 1</p> <p>DigiFocus is a fully digital hearing instrument built around the digital audio processor, DAP. This processor is designed to have a specific functionality and to comply with a series of physical, acoustical and audiological demands. □</p> <p>Oticon – Products Index at Page 1</p> <p> - a 100% digital hearing instrument with fully automatic operation, DigiFocus offers the ultimate in speech understanding. Utilising seven tone band amplification to get the most out of the remaining hearing, DigiFocus is currently the most advanced hearing instrument on the market.</p> <p>Oticon DigiFocus – What Users Say About DigiFocus at p. 1</p> <p>Since January 1996, hundreds of hearing instrument users around the world have participated in field trials of DigiFocus. They all have varying degrees of hearing loss and had worn other instruments for years. During the field trials, Oticon carefully monitored their responses, and the results confirm that DigiFocus is a revolution, not only in terms of technology, but also in improving the quality of life for people who use it.</p> <p>Anderson '783 at Abstract</p> <p>A hearing aid or audio communication system includes an earpiece (10) that can be hidden in the ear canal, and which communicates wirelessly with a remote processor unit, or RPU (16), that enhances audio signals and can be concealed under clothing. Sounds from the environment are picked up by a microphone (12) in the earpiece and sent with other information over a two-way wireless link (17) to the RPU (16). The wireless link (17) uses microwaves for component miniaturization. Furthermore, use of radar technology to implement the wireless link (17), with an RPU (16) interrogator and earpiece (10) transponder, reduces earpiece size and power, as no microwave oscillator is needed in the earpiece (10). Optional secondary wireless link circuitry (19) can be used between the RPU (16) and a cellular telephone system or other sources of information. Electronic voice recognition and response can control system operation.</p>

Claims	Prior Art
	<p data-bbox="583 175 898 207">Anderson '783 at FIG. 1</p>  <p data-bbox="997 747 1060 779">FIG. 1</p>

Claims	Prior Art
	<p data-bbox="583 175 898 207">Anderson '783 at FIG. 2</p>  <p data-bbox="1024 792 1115 824">FIG. 2</p> <p data-bbox="583 846 905 878">Anderson '783 at 1:5–18</p> <p data-bbox="627 883 1948 1133">The present invention is directed to hearing aids, and in particular to a hearing aid having an earpiece housing worn in or at the ear and a remote processor unit (RPU) worn by or located near the user that wirelessly receives signals from and transmits signals to the earpiece. The present invention is also directed to the use of hearing aid systems as wireless communication devices for hands-free cellular telephone and mobile radio communication "handsets," covert operation and control of hearing aids, hearing protection and noise cancellation simultaneous with binaural hearing aid functions, hearing test equipment, location of misplaced hearing aid system components and wireless cochlear implants.</p> <p data-bbox="583 1192 919 1224">Anderson '783 at 1:50–65</p> <p data-bbox="627 1229 1982 1442">The present invention uses an audio path that is different from the audio path used by devices in the prior art. The present invention uses a primary two-way wireless link (not a one-way link as in the prior art) between an earpiece worn in or at the user's ear and an RPU worn by or located near the user. In the present invention, audio signals from the environment are picked up by a microphone in the earpiece (not a hand-held microphone as used by wireless hearing aids in the prior art) and transmitted over the primary two-way wireless link to the RPU (instead of going to processor electronics contained in the earpiece, as</p>

Claims	Prior Art
	<p>in basic hearing aids known in the prior art), where the audio signals are enhanced according to the user's needs before transmission over the primary wireless link to the earpiece. Signal processing is performed in the RPU rather than the earpiece to take advantage of relaxed size and power constraints.</p> <p>Anderson '783 at 3:13–18 FIG. 1 is a block diagram system overview illustrating how an earpiece (comprising a microphone, RF transceiver, dual-use antenna/extractor and speaker) is worn in the ear, communicates via a primary two-way RF link with an RPU, and also communicates via an optional secondary wireless link to a telephone system.</p> <p>Anderson '783 at 3:52–4:14 One preferred embodiment of the invention herein described is shown in block diagram form in FIG. 1. In this figure, the earpiece 10 shown worn in the ear 11 uses a standard completely-in-the-canal (CIC) housing, well known in the art, although many other earpiece housing types known in the art (e.g., behind-the-ear, or BTE) can be used. The earpiece shown in FIG. 1 comprises a microphone 12, RF transceiver 13 with an antenna that doubles as an earpiece extractor 14, and speaker 15. Although the dual-use antenna/extractor arrangement shown in FIG. 1 is often desirable, a separate antenna or antennae can be used in conjunction with many different types of earpiece extractors as appropriate for a particular application. Antennae and extractors can be disguised as ear hair or jewelry (e.g., earrings), or conductive filaments can be permanently implanted in the ear cartilage to act as antennae using a process similar to ear piercing. The earpiece 10 communicates with a remote processor unit (RPU) 16 via a primary two-way RF link 17, although many other wireless link media (e.g., ultrasonic or infrared) may be used instead. The RPU is typically worn under clothing (e.g., carried in a pocket or purse), but may also be worn in plain sight (e.g., on a belt) if desired. The RPU 16 may be connected (via wired or wireless means 18) to optional secondary wireless link circuitry 19 that allows wireless communication between the RPU and other sources of information (e.g., the general subscriber telephone network) via a secondary wireless link. Note that the optional secondary wireless link circuitry 19 may or may not be contained within the RPU case 70 of FIG. 7.</p> <p>Anderson '783 at 5:3–21 During normal operation, speech signals from a nearby talker (and other signals in the ambient audio environment) are picked up by the earpiece wireless microphone 20 and transmitted to the RPU FM receiver 24. The RPU FM receiver 24 output level may be adjusted using the RPU FM receiver 24 volume control. The resulting electrical waveform representing signals from the nearby talker travels through the</p>

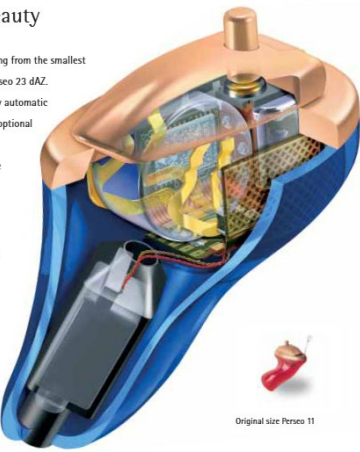
Claims	Prior Art
	<p>DPDT switch 25, which is set to the lower position as shown in FIG. 2, to the signal enhancer 26. The signal enhancer 26 may be, for example, a voice changer device that varies the pitch of a received speech signal according to the setting of pushbutton controls located on the RPU signal enhancer 26. The speech signal's pitch can then be raised or lowered as desired to help compensate for a user's hearing loss relative to a particular talker's voice characteristics. The modified speech signal travels from the RPU signal enhancer 26 to the RPU FM transmitter 27, and finally to the earpiece headset 21 receiver where the signal is converted to acoustic waves heard by the user.</p> <p>Anderson '783 at 22:61–23:35</p> <p>Details of normal system operation: During normal operation of the hearing aid with wireless remote processor, sounds from the ambient environment are picked up by a microphone 12 in the earpiece 10 and transmitted with other information (e.g., data representing the earpiece battery 860 voltage as obtained by an A/D 865 converter and transmitted under control of the earpiece ASIC 845 in the auxiliary bit of many successive replies via the reply data register 825) over a primary two-way wireless link 17 to the RPU 16, where the audio signals are enhanced according to the user's needs by the RPU DSP 948. Signal enhancement is achieved via any of a number of techniques well known in the art that use the capabilities of the general-purpose RPU DSP 948. The enhanced audio signals may be combined with other information (e.g., a synthesized voice, generated within the DSP 948, warning that the aforementioned earpiece battery voltage is low) and transmitted from the RPU 16 over the primary wireless link 17 to the earpiece 10, where they are converted by a speaker 15 to sounds that can only be heard by the user 11. Art auxiliary speaker can be supplied in the RPU 16 as a peripheral 950 device that allows people other than the user to hear the same sounds as those produced by the earpiece speaker 15, if desired (as, e.g., when a person other than the user wishes to use the system's optional cellular telephone capabilities). A different, optional secondary two-way wireless link (shown as part of the RPU in FIG. 9 as secondary wireless link circuitry 944, or connected to the RPU by a wired or other, e.g. infrared, wireless link 18 and shown as optional secondary wireless link circuitry 19 in FIG. 1) can be used for communication between the RPU and a cellular telephone system or other sources of information. An RPU keyboard suitable for covert data entry 72, or voice recognition capabilities in the RPU (implemented, e.g., using any of a number of techniques well known in the art that use the general-purpose RPU DSP 948), can be used to control hearing aid parameters (such as amplification level set by the RPU DSP 948) and telephone dialing functions (e.g., automatic dialing of a predefined telephone number using information stored in the RPU DSP 948) in a manner that can be made imperceptible to a casual observer if desired.</p>

Claims	Prior Art
	<p>Anderson '783 at 27:4–24</p> <p>Hearing test: The present invention allows hearing tests to be performed via the earpiece normally worn by the user, without the need for additional expensive test equipment. All signal generation capabilities necessary to perform hearing tests are available in the hearing aid system comprising an earpiece 10 and RPU 16. Specifically, the RPU DSP 948 can be used to generate tones of varying frequencies and amplitudes, as well as other signals known to be useful for hearing test purposes (e.g., nonsense syllables) that are converted to acoustic waves by the earpiece. Such tones are used to perform a hearing test that determines appropriate gain vs. frequency parameters for a program stored in the RPU DSP 948 that performs signal enhancement to compensate for the user's hearing loss. Note that the hearing test program that controls the RPU DSP 948 during the hearing test can be temporarily stored in the RPU DSP 948 for the duration of the test, then deleted upon completion of the test to allow re-use of RPU DSP 948 memory resources during normal operation. The hearing test program may be loaded into the RPU DSP 948 through the secondary wireless link 944 or a wired peripheral link 950.</p> <p>Enzmann '730 at Abstract</p> <p>An aided ear bud device is provided that operates as a hearing aid device as well as an ear bud device for a wireless phone. The aided ear bud normally operates as a hearing aid device. The aided ear bud state is changed by a verified input from a wireless phone. Upon the change in state, the microphone input from the aided ear bud is remotely output to the wireless phone while the incoming transmission from the wireless phone is combined with the microphone input and output to the user. The ear bud device can be modified to receive analog or digital inputs from other outside sources.</p> <p>Enzmann '730 at 1:63–66</p> <p>Current ear bud devices must be placed in or in close proximity to the ear, meaning that the ear bud placement conflicts with the placement of the hearing aid.</p> <p>Enzmann '730 at 2:17–23</p> <p>Therefore, this invention provides systems and methods for adapting a hearing aid to operate as an ear bud device in addition to providing conventional hearing aid functions. Using this method, this invention is able to address many of the issues hearing impaired wireless phone users encounter while addressing the inherent problems associated with the conventional solutions.</p>

Claims	Prior Art
	<p data-bbox="583 175 915 207">Enzmann '730 at 3:11–14</p> <p data-bbox="632 215 1990 285">The exemplary embodiments of the systems and methods according to this invention provide a method for adapting a hearing aid so that the hearing aid is able to function as an ear bud device i.e., an aided ear bud.</p> <p data-bbox="583 342 915 375">Enzmann '730 at 3:35–40</p> <p data-bbox="632 383 1976 521">The manipulating processes of the DSP 120 serves both to amplify the converted signal from the A/D converter 110 and correct any degraded signal areas so that an improved signal is sent to the speaker 140. The manipulation processes may also adjust the incoming signal based on a preprogrammed frequency mapping of the user's hearing characteristics.</p> <p data-bbox="583 578 915 610">Enzmann '730 at 5:20–30</p> <p data-bbox="632 618 1990 867">The signal 226b passed to the op amp 223 is representative of ambient noise in the aided ear bud 200. This signal 226b is typically reduced in volume and combined with the input from the signal-processing function 225 to produce a mixed signal. The signal processing function 225 processes signal 228 from the TX/RX unit 260 and provides a second processing step so that the timing of the incoming signals from the microphone and the wireless phone are in sync and can be effectively combined to create a true representation of sound to the user. The op amp 223 outputs the mixed signal to the D/A converter 240, and finally to the speaker 250.</p> <p data-bbox="583 924 915 956">Enzmann '730 at 5:34–40</p> <p data-bbox="632 964 1976 1102">The signal 239, output to the speaker 250, is a combined signal based on inputs from the microphone 205 and the wireless phone 270. The combined signal is made to duplicate the experience of a non-hearing impaired person, in that the non-hearing impaired person hears a combination of his or her own voice and ambient noise when speaking into a wireless phone.</p> <p data-bbox="583 1159 1990 1446">To the extent that the DHA systems not disclose any elements of this claim limitation or specific details of implementation, a POSITA would have been motivated to combine elements of DHA system with features known to them and/or disclosed in secondary references, including Anderson or Enzmann, to meet the claim limitations with a reasonable expectation of success. A POSITA would have understood prior art DHA systems to include features and functionality relating to, for example, a personal audio system that generates customized output. It would have been obvious to a POSITA to combine Anderson's and/or Enzmann's teachings for at least the reason that it provides additional information to fill in the missing details of well-known digital hearing aid components and/or functionality. Moreover, such references are readily combinable</p>

Claims	Prior Art
	<p>with DHA systems. For example, the Phonak system, and secondary references like Anderson and Enzmann, share various structural and functional similarities. For example, each discloses ear-worn devices for providing audio output to a user, and include similar components, such as directional microphones, speakers, wireless transceivers, control circuitry, a power source for providing power to the device, and altering audio output. Accordingly, combining the teachings of the Phonak system and Anderson '783 would involve combining known components to yield predictable results.</p> <p>To the extent that the Digital Hearing Aid system and/or Anderson and/or Enzmann does not disclose any elements of this claim limitation, it would have been obvious to a person having ordinary skill in the art to combine, with the teachings of the Digital Hearing Aid system, those elements, as known to them or as disclosed by other prior art identified in the cover pleading.</p>
<p>[1a] a pair of ear buds, each configured to be placed in an ear of a user; and</p>	<p>DHA systems, alone or in view of Anderson, expressly or inherently discloses “a pair of ear buds, each configured to be placed in an ear of a user.”</p> <p>Phonak - SmartLink SX User Guide at p. 32</p> 

Claims	Prior Art
	<p data-bbox="583 175 1010 212">Phonak – Perseo Brochure at p. 8</p> 

Claims	Prior Art
	<p>Phonak – Perseo Brochure at Page 22</p> <h1 data-bbox="655 224 1255 376">The beauty of high technology</h1> <p data-bbox="703 376 961 402">Perseo[®] ITEs – hidden beauty</p> <div data-bbox="630 487 655 535">2</div> <p data-bbox="703 425 1003 587">The Perseo family includes a full line of ITEs ranging from the smallest Perseo 11 CIC to the elegant dual-microphone Perseo 23 dAZ. All CIC and ITE instruments, like the BTEs, are fully automatic and have PersonalLogic. A TacTronic switch is optional on all ITE models and easily added after fitting by exchanging battery covers. This is a very simple and practical means of providing extra control and additional programs.</p> <p data-bbox="703 613 919 695">In many countries, ITEs are available with Phonak eShell. These electronically created shells are hypoallergenic, exceptionally durable and more easily retained in the wearer's ear.</p>  <p data-bbox="1150 787 1234 799">Original size Perseo 11</p>

Claims	Prior Art
	<p data-bbox="583 175 1062 207">Phonak – Perseo Brochure at Page 23</p> <div data-bbox="625 224 1260 824">  <p data-bbox="634 381 903 406">Perseo® BTEs – reliable beauty</p> <p data-bbox="966 438 1260 690">The Perseo BTE family includes a full choice of sizes, colors and power levels. All feature an elegant, award-winning design with a refined, ergonomic shape. They have a patented housing construction with bonded plastic and silicon to effectively keep out debris and moisture. Environmental protection extends further with patented, disposable microphone covers to ensure consistent, reliable performance, in all conditions. With dual-microphone systems this is particularly important since electronic matching cannot successfully compensate for performance reduction as a result of debris and moisture. A comprehensive range of colors for the housing also caters for specific tastes and trends.</p> <p data-bbox="1144 803 1249 820">Original size New small Perseo BTE 111/211</p> </div>

Claims

Prior Art

Phonak – Perseo Brochure at Page 26

Perseo 11

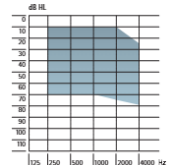


||| Digital microphone

OPTIONS
||| Tactronic switch
||| Control beeps

Ear Simulator data
45 dB max. gain
119 dB SPL max. output

2 cc data
35 dB max. gain
108 dB SPL max. output



Perseo 12

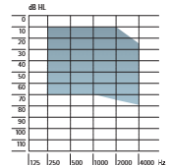


||| Digital microphone
||| Volume control

OPTIONS
||| Tactronic switch
||| Control beeps

Ear Simulator data
50 dB max. gain
119 dB SPL max. output

2 cc data
40 dB max. gain
108 dB SPL max. output



Perseo 22

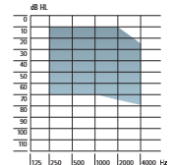


||| Telecoil

OPTIONS
||| Tactronic switch
||| Control beeps
||| Remote control WatchPlot2 and SoundPlot

Ear Simulator data
58 dB max. gain
118 dB SPL max. output

2 cc data
48 dB max. gain
108 dB SPL max. output



Claims

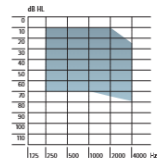
Prior Art

Phonak – Perseo Brochure at Page 27

Perseo 23 dAZ



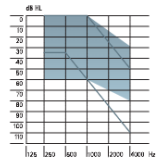
- Adaptive digital AudioZoom
- Telecoil
- OPTIONS
 - Tactronic switch
 - Control beeps
 - Remote control WatchPlot2 and SoundPilot
- Ear Simulator data
 - 58 dB max. gain
 - 118 dB SPL max. output
- 2 cc data
 - 45 dB max. gain
 - 108 dB SPL max. output



Perseo 111 dAZ
Perseo 111 dAZ FM



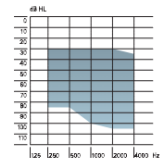
- FM solutions
- Adaptive digital AudioZoom
- Telecoil
- Tactronic switch
- OPTIONS
 - Control beeps
 - Remote control WatchPlot2 and SoundPilot
- Ear Simulator data
 - 54 dB max. gain
 - 127 dB SPL max. output
- 2 cc data
 - 47 dB max. gain
 - 121 dB SPL max. output



Perseo 211 dAZ
Perseo 211 dAZ FM



- FM solutions
- Adaptive digital AudioZoom
- Telecoil
- Tactronic switch
- OPTIONS
 - Control beeps
 - Remote control WatchPlot2 and SoundPilot
- Ear Simulator data
 - 62 dB max. gain
 - 130 dB SPL max. output
- 2 cc data
 - 55 dB max. gain
 - 125 dB SPL max. output



Technical Data

Claims

Prior Art

Phonak – Perseo Brochure at Page 28

Perseo 311 dAZ Forte
Perseo 311 dAZ Forte FM



Remote controls



FM Transmitters



- FM solutions
- Adaptive digital AudioZoom
- Telcoil
- Tachonic Switch

OPTIONS

- Control beeps
- Remote control WatchPilot2 and SoundPilot

Ear Simulator data

72 dB max. gain
135 dB SPL max. output

2 cc data*

85 dB max. gain
130 dB SPL max. output

WatchPilot™2

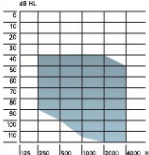


HandyMic™



TelCom™





SoundPilot™



Campus S™



Valente - Effect of type of noise and loudspeaker array on the performance at 398

Differences in performance between omnidirectional and directional microphones were evaluated between two loudspeaker conditions (single loudspeaker at 180°; diffuse using eight loudspeakers set 45° apart) and two types of noise (steady-state HINT noise; R-Space™ restaurant noise). Twenty-five participants were fit bilaterally with Phonak Perseo hearing aids using the manufacturer's recommended procedure. After wearing the hearing aids for one week, the

Claims

Prior Art

Phonak – Collections 2002 Product List at Page 6

Elite


Digital Perception Processing (DPP)


AutoSelect (AS)

New Phonak sound:
unmatched clarity
Digital Perception Processing (DPP) is based on a cochlear sound analysis across 20 critical bands – exactly like the analysis performed by the healthy human ear. But Clari gets a step further than previous overlapping multi-band filter systems. DPP models the active, load-ress-dependent cochlear signal processing in all 20 critical bands.

20 critical bands:
overlapping and dynamic
The cochlear filters are overlapping, but the degree of overlap is dependent on the frequency, level, and bandwidth of the sound signal. This is precisely how Clari's filters work. Digital Perception Processing computes the cochlear excitation patterns and thus determines sound perception. DPP is the first to mimic the healthy human ear in this way.

Automatic, intelligent program selection
Clari listens and automatically activates the optimal DPP algorithm for the acoustic environment: Program 1 with QuietAdapt for quiet situations or Program 2 with NoiseAdapt and noise control for noisy environments.

In quiet situations
QuietAdapt DPP is designed for clear hearing in quiet surroundings. QuietAdapt controls gain in 20 critical bands to correct for loudness and provide a natural sound quality. By preserving the perception patterns of complex sounds, QuietAdapt assures that low-level signal components are not over-

amplified to the point where they would adversely affect speech intelligibility. This allows Clari to deliver only those sound components that would be naturally audible.

In noise
NoiseAdapt DPP is the ideal program for maximizing audibility in noisy environments. The special NoiseAdapt algorithm works with the Fine-scale Noise Canceller – and if available with the Adaptive digital AudioZoom feature – to optimize speech understanding in noise.



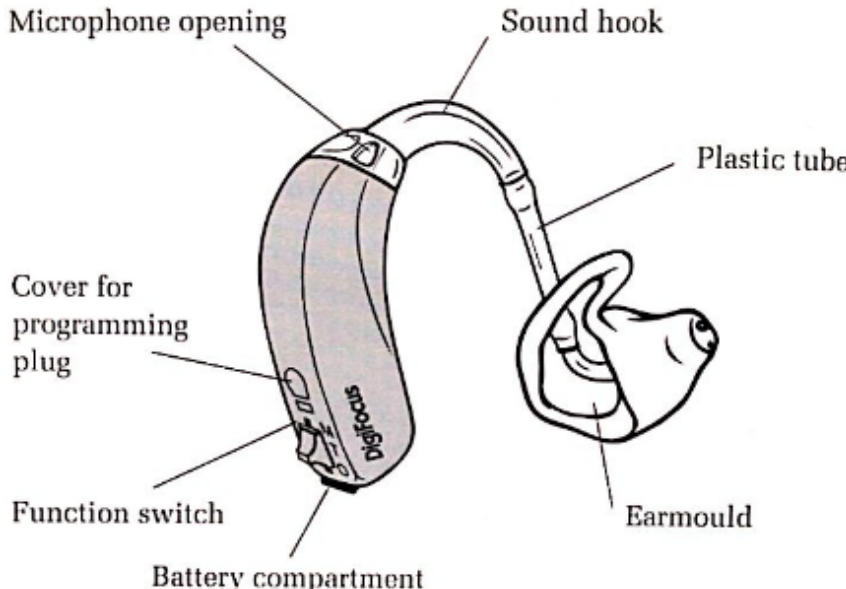
Claims

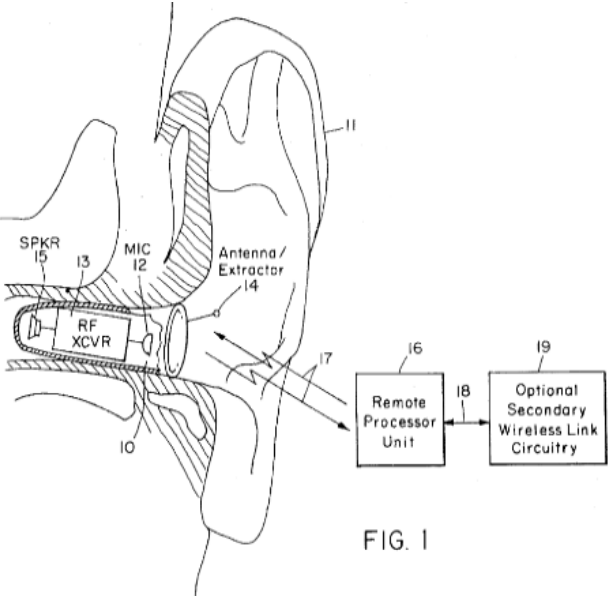
Prior Art

Phonak – Collections 2002 Product List at Page 8

	Claro 11 (CIC)	Claro 21 dAZ	Claro 22
ITE instruments			
			
Peak values G max./OSPL 90	45 dB/117 dB SPL	58 dB/117 dB SPL	58 dB/117 dB SPL
Features	World's first digital microphone TacTronic switch standard	In-the-Ear hearing computer with Adaptive digital Audio/oom (dAZ)	Semi-modular Claro In-the-Canal instrument
Description	Digital Perception Processing (DPP) AutoSelect (AS) Fine-scale Noise Canceler (FNC) 3 hearing programs + "mute" Battery size 10 TacTronic program selection switch	Digital Perception Processing (DPP) AutoSelect (AS) Fine-scale Noise Canceler (FNC) Adaptive digital Audio/oom (dAZ) 3 hearing programs + "mute" Battery size 312 Telecoil	Digital Perception Processing (DPP) AutoSelect (AS) Fine-scale Noise Canceler (FNC) 3 hearing programs + "mute" Battery size 312 Telecoil
Options	Rotary volume control potentiometer Colors: Pink, tan, cocoa, medium brown	Remote controls WatchPlot and SoundPlot FM TacTronic program selection switch On/Off switch (battery cover) Colors: Pink, tan, cocoa, medium brown	Remote controls WatchPlot and SoundPlot FM TacTronic program selection switch On/Off switch (battery cover) Colors: Pink, tan, cocoa, medium brown
Programming	Loudness Perception Profile (LPP) with Active Response Unit (ARU) Fitting formulas: Phonak Digital, Phonak Digital Ski Slope and DSL (i/o)	Loudness Perception Profile (LPP) with Active Response Unit (ARU) Fitting formulas: Phonak Digital, Phonak Digital Ski Slope and DSL (i/o)	Loudness Perception Profile (LPP) with Active Response Unit (ARU) Fitting formulas: Phonak Digital, Phonak Digital Ski Slope and DSL (i/o)

Claims	Prior Art
	<p data-bbox="581 175 1610 212">Lewis – Speech perception in noise: directional microphones (2004) at Page 431</p> <p data-bbox="636 228 1016 818">All subjects were fit with digital Phonak Claro 311 dAZ BTE hearing aids bilaterally. All earmolds had select-a-vent (SAV) venting and #13 or 3 mm horn tubing. In addition to the hearing aids, subjects were fit with Phonak Microlink ML8 FM receivers bilaterally. These FM receivers attach to the bottom of a BTE hearing aid and may be utilized in either the “FM only” mode, which attenuates the hearing aid microphone by 20 dB, or in the FM plus hearing aid mode, which allows for FM input and input of environmental sounds via the hearing aid microphones simultaneously at the same output level. The Phonak TX3 HandyMic FM transmitter served as the FM transmitter. This transmitter has three microphone options: (1) Wide Angle, which picks up sounds arriving from all directions around the transmitter microphone equally; (2) Zoom, which provides reduced amplification to signals arriving from the rear (cardioid); and (3) SuperZoom, which provides reduced amplification to signals originating from the rear and the sides (hypercardioid). The</p>

Claims	Prior Art
	<p data-bbox="583 175 1449 207">Oticon – DigiFocus Compact User Instructional Brochure at Page 4</p>  <p data-bbox="583 844 1081 876">Oticon – DigiFocus Description at p. 1</p> <p data-bbox="630 893 1312 966">DigiFocus is available both as behind-the-ear models (BTE) in colours which match the hair tone, and as in-the-ear (ITE) models which fit discreetly into the ear.</p>  <p data-bbox="714 1291 1239 1315"><i>DigiFocus Compact, Canal, Half Shell and Low Profile.</i></p>

Claims	Prior Art
	<p data-bbox="583 175 1234 207">Oticon – DigiFocus White Paper (1996) at Page 11</p> <p data-bbox="638 232 1310 362">Many people find it uncomfortable to be wired to a computer. This "constricted" feeling can be avoided by using the OtiSet programming belt, which is worn on the shoulder. The programming belt provides cordless infrared communication between the PC and the hearing instrument while in the ear. Alternatively, you may use the HI-PRO interface for a wired solution.</p> <p data-bbox="638 386 1310 540">OtiSet offers a choice of three different views of the fitting. The Audiological View provides a visual guide to instrument performance in four typical sound environments. The Technical View illustrates traditional technical information such as coupler data and insertion gain at different input levels. With the Binaural View you are able to examine settings for both right and left instruments at the same time.</p> <p data-bbox="583 570 894 602">Anderson '783 at FIG. 1</p>  <p data-bbox="583 1235 919 1268">Anderson '783 at 3:13–18</p> <p data-bbox="630 1276 1969 1414">FIG. 1 is a block diagram system overview illustrating how an earpiece (comprising a microphone, RF transceiver, dual-use antenna/extractor and speaker) is worn in the ear, communicates via a primary two-way RF link with an RPU, and also communicates via an optional secondary wireless link to a telephone system.</p>

Claims	Prior Art
	<p data-bbox="583 175 945 207">Anderson '783 at 3:52–4:14</p> <p data-bbox="630 215 1990 833">One preferred embodiment of the invention herein described is shown in block diagram form in FIG. 1. In this figure, the earpiece 10 shown worn in the ear 11 uses a standard completely-in-the-canal (CIC) housing, well known in the art, although many other earpiece housing types known in the art (e.g., behind-the-ear, or BTE) can be used. The earpiece shown in FIG. 1 comprises a microphone 12, RF transceiver 13 with an antenna that doubles as an earpiece extractor 14, and speaker 15. Although the dual-use antenna/extractor arrangement shown in FIG. 1 is often desirable, a separate antenna or antennae can be used in conjunction with many different types of earpiece extractors as appropriate for a particular application. Antennae and extractors can be disguised as ear hair or jewelry (e.g., earrings), or conductive filaments can be permanently implanted in the ear cartilage to act as antennae using a process similar to ear piercing. The earpiece 10 communicates with a remote processor unit (RPU) 16 via a primary two-way RF link 17, although many other wireless link media (e.g., ultrasonic or infrared) may be used instead. The RPU is typically worn under clothing (e.g., carried in a pocket or purse), but may also be worn in plain sight (e.g., on a belt) if desired. The RPU 16 may be connected (via wired or wireless means 18) to optional secondary wireless link circuitry 19 that allows wireless communication between the RPU and other sources of information (e.g., the general subscriber telephone network) via a secondary wireless link. Note that the optional secondary wireless link circuitry 19 may or may not be contained within the RPU case 70 of FIG. 7.</p> <p data-bbox="583 889 919 922">Anderson '783 at 4:30–33</p> <p data-bbox="630 930 1957 995">Although the system described here is monaural, it will be clear to those skilled in the art that a binaural (stereo) system can be constructed based on the same principles as the monaural system.</p> <p data-bbox="583 1052 976 1084">Anderson '783 at 22:61–23:35</p> <p data-bbox="630 1092 1984 1450">Details of normal system operation: During normal operation of the hearing aid with wireless remote processor, sounds from the ambient environment are picked up by a microphone 12 in the earpiece 10 and transmitted with other information (e.g., data representing the earpiece battery 860 voltage as obtained by an A/D 865 converter and transmitted under control of the earpiece ASIC 845 in the auxiliary bit of many successive replies via the reply data register 825) over a primary two-way wireless link 17 to the RPU 16, where the audio signals are enhanced according to the user's needs by the RPU DSP 948. Signal enhancement is achieved via any of a number of techniques well known in the art that use the capabilities of the general-purpose RPU DSP 948. The enhanced audio signals may be combined with other information (e.g., a synthesized voice, generated within the DSP 948, warning that the aforementioned earpiece battery voltage is low) and transmitted from the RPU 16 over the primary wireless link 17 to the</p>

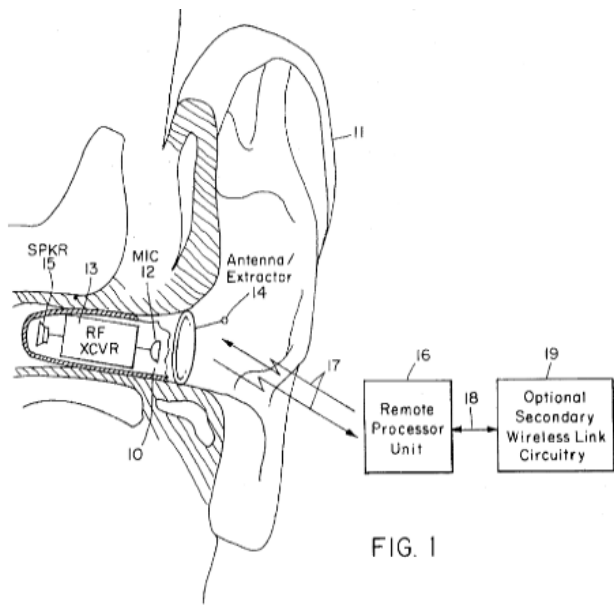
Claims	Prior Art
	<p>earpiece 10, where they are converted by a speaker 15 to sounds that can only be heard by the user 11. Art auxiliary speaker can be supplied in the RPU 16 as a peripheral 950 device that allows people other than the user to hear the same sounds as those produced by the earpiece speaker 15, if desired (as, e.g., when a person other than the user wishes to use the system's optional cellular telephone capabilities). A different, optional secondary two-way wireless link (shown as part of the RPU in FIG. 9 as secondary wireless link circuitry 944, or connected to the RPU by a wired or other, e.g. infrared, wireless link 18 and shown as optional secondary wireless link circuitry 19 in FIG. 1) can be used for communication between the RPU and a cellular telephone system or other sources of information. An RPU keyboard suitable for covert data entry 72, or voice recognition capabilities in the RPU (implemented, e.g., using any of a number of techniques well known in the art that use the general-purpose RPU DSP 948), can be used to control hearing aid parameters (such as amplification level set by the RPU DSP 948) and telephone dialing functions (e.g., automatic dialing of a predefined telephone number using information stored in the RPU DSP 948) in a manner that can be made imperceptible to a casual observer if desired.</p> <p>To the extent that the DHA systems do not disclose any elements of this claim limitation or specific details of implementation, a POSITA would have been motivated to combine elements of the DHA systems with features known to them and/or disclosed in secondary references, including Anderson, to meet the claim limitations with a reasonable expectation of success. A POSITA would have understood prior art DHA systems to include features and functionality relating to, for example, a pair of ear buds, each configured to be placed bilaterally, one in each ear of a user. It would have been obvious to a POSITA to combine Anderson's teachings for at least the reason that it provides information to fill in the missing details of well-known digital hearing aid components and/or functionality, such as speakers, and to improve user experience for users of digital hearing aids seeking to compensate for their hearing impairment. A POSITA would have understood this to also include details relating to directional microphones, wireless transceivers, control circuitry, a power source for providing power to the device, and a mechanical switch for altering audio output. Moreover, such references are readily combinable with DHA systems. For example, the Phonak system, and secondary references like Anderson, share various structural and functional similarities. For example, each discloses ear-worn devices for providing audio output to a user, and include similar components, such as directional microphones, speakers, wireless transceivers, control circuitry, a power source for providing power to the device, and altering audio output. Accordingly, combining the teachings of the Phonak system and Anderson '783 would involve combining known components to yield predictable results.</p> <p>To the extent that the DHA system does not disclose any elements of this claim limitation, it would have been obvious to a person having ordinary skill in the art to combine, with the teachings of DHA systems, those elements, as known to them or as disclosed by other prior art identified in the cover pleading.</p>

Claims	Prior Art
<p>[1b] a mobile communication device at least for supporting audio output, the mobile communication device configured to store at least a hearing profile associated with the user, the stored hearing profile associated with the user at least characterizes hearing impairment regarding the user's ability to hear,</p>	<p>DHA systems, alone or in view of Anderson or Enzmann, expressly or inherently discloses “a mobile communication device at least for supporting audio output, the mobile communication device configured to store at least a hearing profile associated with the user, the stored hearing profile associated with the user at least characterizes hearing impairment regarding the user's ability to hear.”</p> <p>Phonak - SmartLink SX - FM Systems Manual at p. 10</p> <div data-bbox="638 386 1008 418" data-label="Section-Header"> <h4>Smart links with Bluetooth</h4> </div> <div data-bbox="638 423 1092 701" data-label="Text"> <p>With the inclusion of Bluetooth technology in SmartLink, Phonak now enables hearing instrument wearers to enjoy unrestricted use of all Bluetooth compatible cellular phones. One of the major complaints of hearing impaired people is now history – talking on the telephone is even easier for them than for people with normal hearing. Simply synchronizing the Bluetooth telephone and SmartLink creates the communication path. As long as the cell phone is within a 10-meter range of SmartLink, all incoming calls are automatically directed</p> </div> <div data-bbox="1144 396 1600 701" data-label="Text"> <p>via SmartLink to the hearing instrument. The cell phone is only needed for dialing outgoing calls. Whenever the cell phone is switched on, the Bluetooth function remains active for incoming calls, regardless of whether SmartLink is switched on or off. This guarantees that telephone calls are not missed, even when the FM function is not required. When using an external audio input to listen to music, incoming telephone calls have priority when the cell phone is switched on. The audio signal is automatically decreased so that the user hears the ringing tone clearly.</p> </div> <p>Phonak – Perseo Brochure at Page 5</p> <div data-bbox="638 781 1247 834" data-label="Section-Header"> <h3>The Science of PersonalLogic™</h3> </div> <div data-bbox="741 886 1793 1000" data-label="Text"> <p>Phonak has applied knowledge from new research and past experience to add a new dimension of individualized logic to automatic hearing systems – PersonalLogic. Combining new and proven features, PersonalLogic optimizes the functioning, operation and fitting of hearing systems for each individual user.</p> </div> <div data-bbox="741 1058 1820 1258" data-label="Text"> <p>Thanks to PersonalLogic, Perseo takes customization to a level beyond “custom-fitted” related to individual hearing loss. The new digital system from Phonak is the first hearing system with customized automatic operation, as opposed to operation based on factory set, unchangeable rules. Now, Perseo gives hearing professionals control over automatic operation, allowing them to provide clients with a hearing system “custom-tailored” to their personal hearing requirements and tastes.</p> </div>

Claims	Prior Art
	<p>Phonak – Perseo Brochure at Page 7</p> <div data-bbox="636 232 737 332" style="display: inline-block; width: 40px; height: 60px; background-color: #FFD700; margin-right: 10px;"></div> <div data-bbox="758 228 1207 267" style="display: inline-block;">Convenience and Comfort</div> <div data-bbox="758 271 1100 295" style="display: inline-block;">Perseo is easier to wear and easier to fit</div> <div data-bbox="732 354 1224 636" style="display: inline-block; vertical-align: top;"> <ul style="list-style-type: none"> ■ For true "wear and hear" convenience and comfort in all situations, Perseo automatic operation is managed by Phonak-exclusive PersonalLogic. ■ This breakthrough intelligence for automatic hearing instruments responds to acoustics, the listening situation and personal listening preferences. ■ Perseo is easily programmed for individual needs using the Listening Situation Manager. </div> <div data-bbox="1230 232 1331 332" style="display: inline-block; width: 40px; height: 60px; background-color: #00AEEF; margin-right: 10px;"></div> <div data-bbox="1352 228 1470 267" style="display: inline-block;">Clarity</div> <div data-bbox="1352 271 1816 295" style="display: inline-block;">Enhanced performance in multiple listening situations</div> <div data-bbox="1329 354 1820 750" style="display: inline-block; vertical-align: top;"> <ul style="list-style-type: none"> ■ The Clarity Component Enhancement system realized in Digital Perception Processing² (DPP²) increases the clarity of speech and environmental sounds. ■ Rapid-tracking Adaptive digital AudioZoom (dAZ) is automatically activated to improve speech understanding in noisy situations. ■ For extra comfort and ease of listening, the Fine-scale Noise Canceler selectively reduces background noise in 20 channels. ■ Optional MicroLink FM provides the ultimate Speech-to-Noise Ratio enhancement. </div> <p>Phonak – Perseo Brochure at p. 17</p> <p>With the development of Perseo, PFG has become significantly easier to use, providing direct access to the benefits of high-tech hearing systems. PFG for Perseo pre-sets the instrument according to each wearer's hearing ability and needs. The experience gained over a decade of PFG fittings and extensive field trials is integrated automatically into initial settings. This ensures maximum spontaneous acceptance.</p> <p>The wearer immediately feels confident that Perseo is the right choice. This confidence continues even if additional adjustments are necessary, thanks to the flexibility of the three levels of professional control – Easy with Personal System Managers, Standard with clustered controls, and Detailed with access to each of the 20 channels.</p>

Claims	Prior Art
	<p>Valente - Effect of type of noise and loudspeaker array on the performance at 402</p> <p>The Phonak fitting software (Version 8.3) was used to program the hearing aids. First, hearing thresholds were entered into the NOAH software. The manufacturer's prescribed "First-Fit" was calculated from these hearing thresholds. The default settings in each of the software "personal system" menus (provides adjustments to the "First-Fit" for prior experience; occlusion effect and feedback) were utilized. If the patient reported occlusion, then this feature was activated. If feedback was noticed, the feedback manager was accessed. All hearing</p> <p>Phonak – Online Interview re Perseo (April 7, 2003) at Page 7</p> <p>PersonalLogic may be personalized, via the Listening Situation Manager, to meet individual patient preferences for listening in noise. The clinician may select whether PersonalLogic adapts to noise on the basis of audibility or comfort. Audibility changes to Noisy situation processing only when speech in noise is detected, but not with speech alone or noise alone. Comfort adapts WHENEVER noise is present that exceeds a certain dB level. Are you with me?</p> <p>Phonak – Fitting Guide for MLxi and MLxi Baha at Page 8</p> <div data-bbox="630 901 1606 1453"> <div data-bbox="630 901 1113 1453"> <p>7. Programming options</p> <p>7.1 Programming interface</p>  <p>All receivers can be programmed using the FM programming interface:</p> <ol style="list-style-type: none"> 1. Make sure all equipment is charged or a fresh battery is inserted in the hearing instrument. Turn all equipment on. 2. Connect the universal receiver to the hearing instrument and insert it into the programming interface. 3. Always select the option "via the programming interface" when FM SuccessWare requests communication with the receiver. </div> <div data-bbox="1113 901 1606 1453"> <p>7.2 Using inspiro</p>  <p>MLxi, MLxi Baha can be programmed via inspiro:</p> <ol style="list-style-type: none"> 1. Make sure a fresh battery is inserted in the hearing instrument. Turn all equipment on. 2. Connect the universal receiver to the hearing instrument, the CI speech processor, or MicroLink CI S₂, and place in front of the inspiro. 3. Always select the option "via inspiro" when the FM SuccessWare requests communication with the receiver. </div> </div>

Claims	Prior Art
	<p data-bbox="583 175 1234 207">Phonak – Collections 2002 Product List at Page 33</p> <div data-bbox="636 256 804 337"> <p>PFG – new perspectives in fitting success</p> </div> <div data-bbox="636 386 884 570">  </div> <div data-bbox="636 591 884 773">  </div> <div data-bbox="945 228 1297 282"> <p>QuickGuide The versatile, easy-to-use online tour through the PFG fitting process with access to all optional fitting tools.</p> </div> <div data-bbox="945 303 1297 355"> <p>Instrument selection wizard This online assistant recommends the product that best matches the customer's needs.</p> </div> <div data-bbox="945 376 1297 449"> <p>Initial SoundCheck This is an ideal tool for fast and simple fine tuning of sound quality. Three controls provide convenient gain and balance settings in the high- and low-frequency ranges.</p> </div> <div data-bbox="945 470 1297 524"> <p>Automatic fine tuning Based on the customer's hearing perception, the software automatically computes the most suitable modifications.</p> </div> <div data-bbox="945 545 1297 599"> <p>Manual fine tuning A high degree of fitting precision can be attained with this broad spectrum of parameters in an optimised display format.</p> </div> <div data-bbox="945 620 1297 693"> <p>Phonak Sound The degree of fitting success can be dramatically improved by using these real-life sound clips. All sound examples can be saved on a computer as MP3 files; they are also available on a sound CD.</p> </div> <div data-bbox="945 714 1197 820"> <p>Proven PFG benefits</p> <ul style="list-style-type: none"> • One program for all Phonak hearing computers • NOAH compatible or standalone version • Detailed help system with in-depth scientific and technical information • PFG training demo </div>

Claims	Prior Art
	<p data-bbox="583 175 898 207">Anderson '783 at FIG. 1</p>  <p data-bbox="997 747 1060 779">FIG. 1</p>

Claims	Prior Art
	<p>in basic hearing aids known in the prior art), where the audio signals are enhanced according to the user's needs before transmission over the primary wireless link to the earpiece. Signal processing is performed in the RPU rather than the earpiece to take advantage of relaxed size and power constraints.</p> <p>Anderson '783 at 2:60–3:5 A further object of this invention is to provide some degree of protection for the residual hearing capability of a hearing-impaired user (or a non-impaired user who wishes to avoid impairment) in a wireless hearing aid system having the simultaneous capabilities of noise cancellation and binaural processing (e.g., directionality). Other objects of the invention are to provide a convenient means for testing a user's hearing capability without the need for additional equipment, to assist users in the location of misplaced hearing aid system components, and to provide profoundly deaf cochlear implant patients with a wireless system allowing improved freedom of movement compared to existing wire-connected systems.</p> <p>Anderson '783 at 3:42–43 FIG. 9 is a block diagram showing details of an RPU for a full-featured preferred embodiment of the invention.</p> <p>Anderson '783 at 5:3–21 During normal operation, speech signals from a nearby talker (and other signals in the ambient audio environment) are picked up by the earpiece wireless microphone 20 and transmitted to the RPU FM receiver 24. The RPU FM receiver 24 output level may be adjusted using the RPU FM receiver 24 volume control. The resulting electrical waveform representing signals from the nearby talker travels through the DPDT switch 25, which is set to the lower position as shown in FIG. 2, to the signal enhancer 26. The signal enhancer 26 may be, for example, a voice changer device that varies the pitch of a received speech signal according to the setting of pushbutton controls located on the RPU signal enhancer 26. The speech signal's pitch can then be raised or lowered as desired to help compensate for a user's hearing loss relative to a particular talker's voice characteristics. The modified speech signal travels from the RPU signal enhancer 26 to the RPU FM transmitter 27, and finally to the earpiece headset 21 receiver where the signal is converted to acoustic waves heard by the user.</p> <p>Anderson '783 at 22:61–23:35 Details of normal system operation: During normal operation of the hearing aid with wireless remote</p>

Claims	Prior Art
	<p>processor, sounds from the ambient environment are picked up by a microphone 12 in the earpiece 10 and transmitted with other information (e.g., data representing the earpiece battery 860 voltage as obtained by an A/D 865 converter and transmitted under control of the earpiece ASIC 845 in the auxiliary bit of many successive replies via the reply data register 825) over a primary two-way wireless link 17 to the RPU 16, where the audio signals are enhanced according to the user's needs by the RPU DSP 948. Signal enhancement is achieved via any of a number of techniques well known in the art that use the capabilities of the general-purpose RPU DSP 948. The enhanced audio signals may be combined with other information (e.g., a synthesized voice, generated within the DSP 948, warning that the aforementioned earpiece battery voltage is low) and transmitted from the RPU 16 over the primary wireless link 17 to the earpiece 10, where they are converted by a speaker 15 to sounds that can only be heard by the user 11. Art auxiliary speaker can be supplied in the RPU 16 as a peripheral 950 device that allows people other than the user to hear the same sounds as those produced by the earpiece speaker 15, if desired (as, e.g., when a person other than the user wishes to use the system's optional cellular telephone capabilities). A different, optional secondary two-way wireless link (shown as part of the RPU in FIG. 9 as secondary wireless link circuitry 944, or connected to the RPU by a wired or other, e.g. infrared, wireless link 18 and shown as optional secondary wireless link circuitry 19 in FIG. 1) can be used for communication between the RPU and a cellular telephone system or other sources of information. An RPU keyboard suitable for covert data entry 72, or voice recognition capabilities in the RPU (implemented, e.g., using any of a number of techniques well known in the art that use the general-purpose RPU DSP 948), can be used to control hearing aid parameters (such as amplification level set by the RPU DSP 948) and telephone dialing functions (e.g., automatic dialing of a predefined telephone number using information stored in the RPU DSP 948) in a manner that can be made imperceptible to a casual observer if desired.</p> <p>Anderson '783 at 27:4–24</p> <p>Hearing test: The present invention allows hearing tests to be performed via the earpiece normally worn by the user, without the need for additional expensive test equipment. All signal generation capabilities necessary to perform hearing tests are available in the hearing aid system comprising an earpiece 10 and RPU 16. Specifically, the RPU DSP 948 can be used to generate tones of varying frequencies and amplitudes, as well as other signals known to be useful for hearing test purposes (e.g., nonsense syllables) that are converted to acoustic waves by the earpiece. Such tones are used to perform a hearing test that determines appropriate gain vs. frequency parameters for a program stored in the RPU DSP 948 that performs signal enhancement to compensate for the user's hearing loss. Note that the hearing test program that controls the RPU DSP 948 during the hearing test can be temporarily stored in the RPU DSP 948 for the duration of the test, then deleted upon completion of the test to allow re-use of RPU DSP 948 memory resources during normal operation. The hearing test program may be loaded into the RPU DSP 948</p>

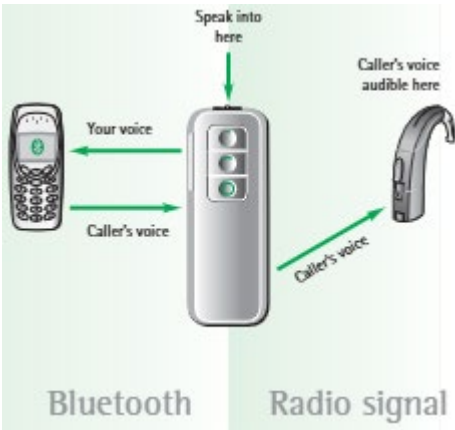
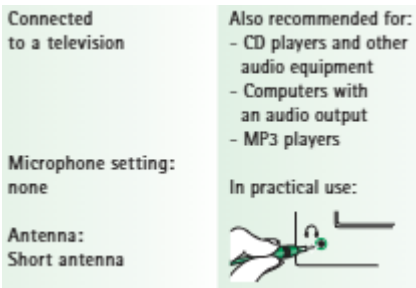
Claims	Prior Art
	<p>through the secondary wireless link 944 or a wired peripheral link 950.</p> <p>Similarly, Enzmann '730 discloses storing at least a hearing profile associated with the user that characterizes a user's hearing impairment.</p> <p>Enzmann '730 at 6:1–11</p> <p>The digital logic 320 contains a series of devices consisting of multiple identical and separately identifiable structures built of some mix of combinational logic and sequential elements the sequential elements being commonly referred to as registers, latches, flip-flops, or shift register latches. These devices individually or in combination implement a series of activation states and wave patterns that determine the activation state of the aided ear bud and shapes of the frequency wave output from the digital logic 320. The functionality of the digital logic 320 is depicted in FIG. 4 as the amplification chart 325.</p> <p>Enzmann '730 at 6:12–35</p> <p>The amplification chart 325 contains a series of pre-stored dynamic frequency maps based on the user's hearing characteristics, the speaker characteristics and the ambient noise near the aided ear bud 200. The digital logic 320 will receive the transmitted signal from the TX/RX unit 260 and process those signals to an optimal form based on the user's hearing, speaker and ambient noise characteristics. The arrangement of the logic devices within the digital logic 320 will serve to produce each of the characteristics (wave shape) of the amplitude chart 325. The user's hearing characteristics remain static unless the aided ear bud 200 is reprogrammed to compensate for additional hearing loss or increase and different users. For example, as shown in FIG. 5, a mapping of a user's hearing characteristic (denoted by the “••” line) illustrates that the user lacks significant hearing capability in higher frequencies. The desired hearing characteristic is denoted by the “-••” line. This desired hearing characteristic depicts or represents how a non-hearing impaired person's hearing characteristics may be. As known to those skilled in the art, close representation of this hearing characteristic is achievable by adding a complementary characteristic to the user's hearing characteristic to produce a result that is very close to the desired characteristic. Each of the characteristics is corrected using the process shown in FIG. 5.</p>

Claims	Prior Art
	<p data-bbox="583 175 894 207">Enzmann '730 at FIG. 5</p> <p data-bbox="953 800 1033 833">FIG. 5</p> <p data-bbox="583 849 915 881">Enzmann '730 at 6:44–48</p> <p data-bbox="632 886 1992 992">It should be appreciated that the characteristics represented in FIG. 5 are for illustrative purposes only. Each of the characteristics in FIG. 5 are device and user specific. Each of the characteristics can be set and calculated by information stored in the EEPROM 310.</p> <p data-bbox="583 1011 1992 1446">A POSITA would have been motivated to combine the DHA system with the features of Anderson and/or Enzmann to meet the claim limitations with a reasonable expectation of success. For example, a POSITA would have understood the disclosures of the Phonak system to teach well-known structures and functionality associated with a digital hearing aid device, including storing a hearing profile/characteristics for the user based on the results of the hearing test, i.e., such that the hearing profile/characteristics stores the appropriate hearing parameters for the user on the ear-level device itself, and using the hearing profile/characteristics to compensate for the user's hearing loss, as described above in the citations. A POSITA would have found it obvious to configure and use an external device, such a mobile communication device, with the functionality necessary to store a user's hearing profile, execute algorithms for customizing audio according to that stored hearing profile, and transmit the processed audio back to the hearing aid device, in order to leverage the external device's existing communication hardware and processing capabilities to expand the digital hearing aid's functionality to audio signals generated or received by the external device itself, before or instead of</p>

Claims	Prior Art
	<p>signals received by the digital hearing aid itself. To the extent the specific details of features/functionality of similar digital hearing aids are not expressly disclosed, a POSITA thus would have been motivated to look to other references, like Anderson or Enzmann, to supply these missing details. A POSITA would have understood that enabling the mobile communication device with such processing capabilities would reduce power burden and antenna requirements of the hearing aid device and provide an improved user experience by reducing signal drop outs that users frequently experience with signal being transmitted from the external device for processing by the digital hearing aid. A POSITA would have understood that this approach would eliminate the need for large and power-consuming electronics in the earpiece. A POSITA would have had a reasonable expectation of success in implementing this combination, as wireless communication between host and ear-level devices and digital signal processing (on the ear-level device for customization according to a user's hearing profile and on an external host device for enhancement and acoustic feedback control) were well-established practices at the relevant time. The combination therefore would have required no more than the routine application of known engineering techniques and would have predictably achieved the intended result to configuring the remote device with functionality already existing on the digital hearing aid.</p> <p>To the extent that the DHA systems and Anderson or Enzmann do not disclose any elements of this claim limitation, it would have been obvious to a person having ordinary skill in the art to combine, with the teachings of DHA systems, those elements, as known to them or as disclosed by other prior art identified in the cover pleading.</p>

Claims	Prior Art
<p>[1c] the mobile communication device configured to produce audio output signals customized for the user based on at least the stored hearing profile associated with the user such that the customized audio output signals at least compensate for the hearing impairment regarding the user's ability to hear,</p>	<p>DHA systems in view of Anderson, expressly or inherently discloses “the mobile communication device configured to produce audio output signals customized for the user based on at least the stored hearing profile associated with the user such that the customized audio output signals at least compensate for the hearing impairment regarding the user's ability to hear.”</p> <p>Phonak - SmartLink SX - FM Systems Manual at p. 3</p> <p>SmartLink, based on Phonak breakthrough digital FM synthesizer technology, now enables people with a hearing impairment and those with normal hearing to use the same communication tools with equal effectiveness. SmartLink is a communication gateway to all situations where understanding perfectly, all of the time, is key. The unique innovation, SmartLink, meets the needs of people with hearing loss, giving them access to effective communication in difficult situations and to the many advantages off modern communication tools, via one device.</p> <p>Phonak - SmartLink SX - FM Systems Manual at p. 4</p> <p>Smart connections to all audio sources Any external audio input or an external microphone signal can be transmitted wirelessly via SmartLink to the hearing instrument, simply by attaching to the smart interface.</p> <p>Phonak - SmartLink SX - FM Systems Manual at p. 10</p> <p>Smart links with Bluetooth With the inclusion of Bluetooth technology in SmartLink, Phonak now enables hearing instrument wearers to enjoy unrestricted use of all Bluetooth compatible cellular phones. One of the major complaints of hearing impaired people is now history – talking on the telephone is even easier for them than for people with normal hearing. Simply synchronizing the Bluetooth telephone and SmartLink creates the communication path. As long as the cell phone is within a 10-meter range of SmartLink, all incoming calls are automatically directed via SmartLink to the hearing instrument. The cell phone is only needed for dialing outgoing calls. Whenever the cell phone is switched on, the Bluetooth function remains active for incoming calls, regardless of whether SmartLink is switched on or off. This guarantees that telephone calls are not missed, even when the FM function is not required. When using an external audio input to listen to music, incoming telephone calls have priority when the cell phone is switched on. The audio signal is automatically decreased so that the user hears the ringing tone clearly.</p>

Claims	Prior Art
	<p data-bbox="583 175 1289 211">Phonak - SmartLink SX - FM Systems Manual at p. 10</p>  <p data-bbox="583 841 1150 876">Phonak - SmartLink SX User Guide at p. 27</p> <ul data-bbox="630 893 1050 974" style="list-style-type: none"> • The SmartLink can be directly connected to audio equipment like televisions, CD and MP3 players, and computers with audio outputs.

Claims	Prior Art
	<p data-bbox="583 175 1150 211">Phonak - SmartLink SX User Guide at p. 32</p>  <p data-bbox="583 670 1150 706">Phonak - SmartLink SX User Guide at p. 47</p> 

Claims

Prior Art

Phonak – Perseo Brochure at Page 15

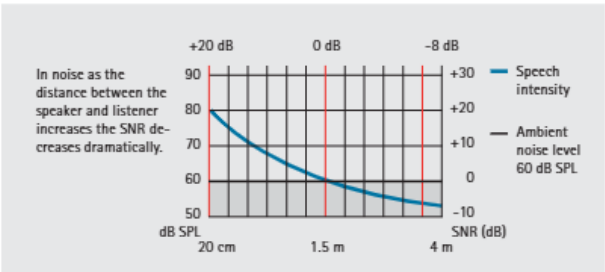
The perfectly integrated communication system

FM MicroLink

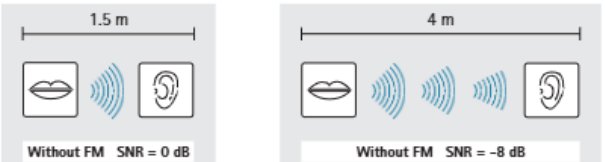
FM technology has allowed thousands of children to participate in standard education settings and to be actively involved in social activities. Phonak pioneered this development and the integration of FM into hearing instruments to produce communication systems that are effective in even the most demanding situations. Since the innovation of miniaturized Phonak **MicroLink** systems, FM has become an important solution also for adults with a high communication need in work and leisure environments.

All Perseo BTEs have been designed to incorporate MicroLink with its fully integrated appearance and operation. This integrated system is unique to Phonak.

Speech-to-Noise Ratio (SNR) as a function of distance between speaker and listener



Speech-to-Noise Ratio without FM is distance dependent

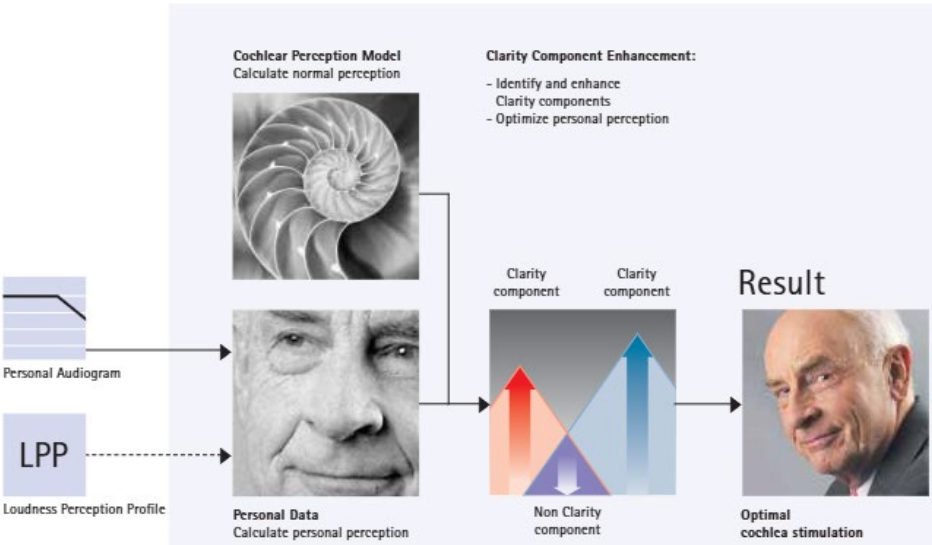


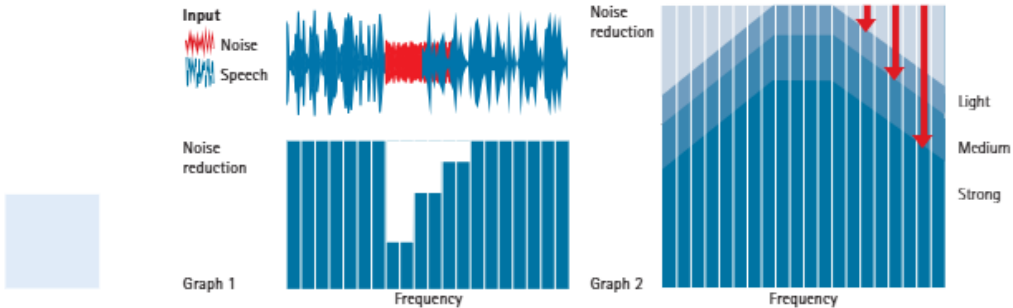
Claims	Prior Art
	<p data-bbox="583 175 1224 207">Phonak SmartLink SX Product Information at p. 1</p> <div data-bbox="630 224 1125 816"> <p data-bbox="642 228 856 253">Product Information</p> <p data-bbox="642 264 980 315">SmartLink™ SX</p>  <p data-bbox="642 659 1121 816">SmartLink SX is a compact communication device for people with hearing loss who lead busy, active lives. The SmartLink SX incorporates an FM transmitter with digital signal processing, a remote control for Phonak hearing instruments, and Bluetooth for communication with mobile phones. With its digital frequency synthesizer, the SmartLink allows a highly flexible choice of operating frequency. For maximum customer benefit, we recommend using a Phonak multi-channel FM receiver (e.g. the MLxS or MLBS) that can be automatically synchronized from the SmartLink.</p> </div> <p data-bbox="583 846 1491 878">Phonak – Media Release re Product Offerings in F2003/2004 at Page 2</p> <p data-bbox="636 894 1518 1114">The key breakthrough technology in Perseo is <i>PersonalLogic</i>, the first ever automatic control system that is fitter programmable according to the tastes and preferences of each individual wearer. Automatic performance has now a new level of intelligence beyond merely adapting to acoustic changes in the environment. Perseo is the first hearing system with customized automatic operation, as opposed to operation based on factory settings. It allows hearing professionals to provide clients with a hearing system "custom-tailored" to their personal hearing requirements and tastes.</p> <p data-bbox="583 1143 974 1175">Anderson '783 at 22:61–23:35</p> <p data-bbox="630 1179 1982 1468">Details of normal system operation: During normal operation of the hearing aid with wireless remote processor, sounds from the ambient environment are picked up by a microphone 12 in the earpiece 10 and transmitted with other information (e.g., data representing the earpiece battery 860 voltage as obtained by an A/D 865 converter and transmitted under control of the earpiece ASIC 845 in the auxiliary bit of many successive replies via the reply data register 825) over a primary two-way wireless link 17 to the RPU 16, where the audio signals are enhanced according to the user's needs by the RPU DSP 948. Signal enhancement is achieved via any of a number of techniques well known in the art that use the capabilities of the general-purpose RPU DSP 948. The enhanced audio signals may be combined with other</p>

Claims	Prior Art
	<p>information (e.g., a synthesized voice, generated within the DSP 948, warning that the aforementioned earpiece battery voltage is low) and transmitted from the RPU 16 over the primary wireless link 17 to the earpiece 10, where they are converted by a speaker 15 to sounds that can only be heard by the user 11. Art auxiliary speaker can be supplied in the RPU 16 as a peripheral 950 device that allows people other than the user to hear the same sounds as those produced by the earpiece speaker 15, if desired (as, e.g., when a person other than the user wishes to use the system's optional cellular telephone capabilities). A different, optional secondary two-way wireless link (shown as part of the RPU in FIG. 9 as secondary wireless link circuitry 944, or connected to the RPU by a wired or other, e.g. infrared, wireless link 18 and shown as optional secondary wireless link circuitry 19 in FIG. 1) can be used for communication between the RPU and a cellular telephone system or other sources of information. An RPU keyboard suitable for covert data entry 72, or voice recognition capabilities in the RPU (implemented, e.g., using any of a number of techniques well known in the art that use the general-purpose RPU DSP 948), can be used to control hearing aid parameters (such as amplification level set by the RPU DSP 948) and telephone dialing functions (e.g., automatic dialing of a predefined telephone number using information stored in the RPU DSP 948) in a manner that can be made imperceptible to a casual observer if desired.</p> <p>Anderson '783 at 27:4–24</p> <p>Hearing test: The present invention allows hearing tests to be performed via the earpiece normally worn by the user, without the need for additional expensive test equipment. All signal generation capabilities necessary to perform hearing tests are available in the hearing aid system comprising an earpiece 10 and RPU 16. Specifically, the RPU DSP 948 can be used to generate tones of varying frequencies and amplitudes, as well as other signals known to be useful for hearing test purposes (e.g., nonsense syllables) that are converted to acoustic waves by the earpiece. Such tones are used to perform a hearing test that determines appropriate gain vs. frequency parameters for a program stored in the RPU DSP 948 that performs signal enhancement to compensate for the user's hearing loss. Note that the hearing test program that controls the RPU DSP 948 during the hearing test can be temporarily stored in the RPU DSP 948 for the duration of the test, then deleted upon completion of the test to allow re-use of RPU DSP 948 memory resources during normal operation. The hearing test program may be loaded into the RPU DSP 948 through the secondary wireless link 944 or a wired peripheral link 950.</p> <p>A POSITA would have been motivated to combine the teachings of the DHA system with Anderson to meet the claim limitations with a reasonable expectation of success for the reasons described above with respect to claim limitation [1b]. A POSITA would have been further motivated to combine the DHA system with the features of Anderson in order to leverage the external device's existing communication hardware and</p>

Claims	Prior Art
	<p>processing capabilities to expand the digital hearing aid’s functionality to audio signals generated or received by the external device itself, before or instead of signals received by the digital hearing aid itself. A POSITA would have understood the disclosures of the Phonak or other DHA system to teach well-known structures and functionality associated with a digital hearing aid device, including storing a hearing profile/characteristics for the user based on the results of the hearing test, i.e., such that the hearing profile/characteristics stores the appropriate hearing parameters for the user on the ear-level device itself, and using the hearing profile/characteristics to compensate for the user’s hearing loss, as described above in the citations. Configuring and using an external device, such a mobile communication device, with the functionality necessary to store a user’s hearing profile, execute algorithms for customizing audio according to that stored hearing profile, and transmit the processed audio back to the hearing aid device, would allow those with hearing impairments to continue using their DHA system simultaneously with their mobile communication device without sacrificing the benefits provided by the DHA system. A POSITA would have also understood that enabling the mobile communication device with such processing capabilities would the reduce power burden and antenna requirements of the hearing aid device and provide an improved user experience by reducing signal drop outs that users frequently experience with signal being transmitted from the external device for processing by the digital hearing air, and eliminate the need for large and power-consuming electronics in the earpiece. A POSITA would have had a reasonable expectation of success in implementing this combination, as wireless communication between host and ear-level devices and digital signal processing (on the ear-level device for customization according to a user’s hearing profile and on an external host device for enhancement and acoustic feedback control) were well-established practices at the relevant time. The combination therefore would have required no more than the routine application of known engineering techniques and would have predictably achieved the intended result to configuring the remote device with functionality already existing on the digital hearing aid.</p> <p>To the extent that the DHA system and Anderson not disclose any elements of this claim limitation, it would have been obvious to a person having ordinary skill in the art to combine, with the teachings of the DHA system, those elements, as known to them or as disclosed by other prior art identified in the cover pleading.</p>
<p>[1d] the mobile communication device configured to produce the customized audio output signals using at least signal amplification, and</p>	<p>DHA systems, alone or in view of Anderson, expressly or inherently discloses “the mobile communication device configured to produce the customized audio output signals using at least signal amplification.”</p>

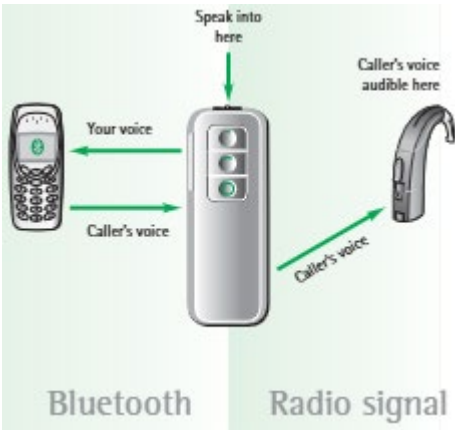
Claims	Prior Art
	<p>Phonak - SmartLink SX - FM Systems Manual at p. 4</p> <div data-bbox="632 228 1098 470"> <p>Smart signal processing For the first time in the industry, SmartLink combines all of the verified advantages of digital signal processing strategies with the proven Phonak miniaturized MicroLink technology.</p> <p>SmartLink – the gateway to crystal-clear sound quality</p> </div> <p>Phonak - SmartLink SX - FM Systems Manual at p. 6</p> <div data-bbox="632 540 917 688"> <p>Advanced DSP processing package:</p> <ul style="list-style-type: none"> • Digital SuperZoom² • Fine-scale Noise Canceler • Digital voice enhancement </div> <p>Phonak - SmartLink SX - FM Systems Manual at p. 9</p> <p>The SmartLink three-phase signal recognition and audio processing strategy is the key to creating natural, crystal-clear sound. In the first phase, the focus of reception is defined through the choice of microphone settings. Secondly, the noise cancellation eliminates disturbing noises detected within the defined beam. Finally, voices within the defined beam are enhanced through digital voice enhancement technology.</p> <p>Phonak - SmartLink SX - FM Systems Manual at p. 9</p> <p>Digital voice enhancement By analyzing the acoustic characteristics of the incoming sound, the digital voice enhancement feature optimizes the perceived loudness of the transmitted FM signal.</p>

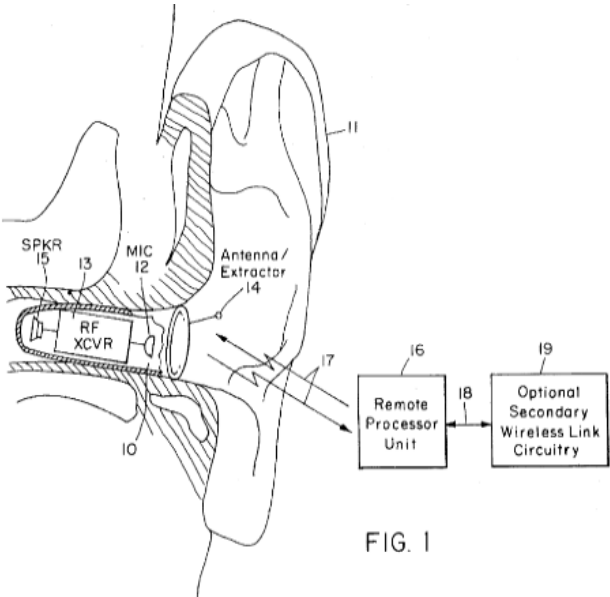
Claims	Prior Art
	<p data-bbox="583 175 1062 207">Phonak – Perseo Brochure at Page 13</p> <div data-bbox="625 264 1482 310"> <h2>Enhanced sound quality and Speech-to-Noise Ratio</h2>  </div> <p data-bbox="726 362 1062 383">Clarity Component Enhancement – DPP²</p> <p data-bbox="726 401 1719 617">Speech clarity depends on both its acoustic properties and the listener's auditory perception. Conventional spectral enhancement systems consider only the acoustics of sound, not how it will be perceived by the listener. Perseo goes well beyond this with DPP², the new generation of Digital Perception Processing. DPP² with Clarity Component Enhancement is individualized for each wearer. It analyzes the acoustic and the psychoacoustic characteristics of incoming sound in 20 critical bands to identify accurately the components that contribute to clarity. It then selectively amplifies them to optimize clarity for the wearer.</p> <p data-bbox="583 638 1062 670">Phonak – Perseo Brochure at Page 13</p> <div data-bbox="636 686 1169 725"> <h3>Digital Perception Processing² (DPP²)</h3>  <pre> graph LR PA[Personal Audiogram] --> CPM[Cochlear Perception Model Calculate normal perception] LPP[LPP Loudness Perception Profile] -.-> PD[Personal Data Calculate personal perception] CPM --> CCE[Clarity Component Enhancement: - Identify and enhance Clarity components - Optimize personal perception] PD --> CCE CCE --> R[Result Optimal cochlea stimulation] CCE --> NC[Non Clarity component] CCE --> CC1[Clarity component] CCE --> CC2[Clarity component] </pre> </div>

Claims	Prior Art
	<p data-bbox="583 175 1064 207">Phonak – Perseo Brochure at Page 14</p> <p data-bbox="632 228 1205 261">Enhanced comfort in all situations</p> <p data-bbox="722 272 1073 297">Fine-scale Noise Canceler in 20 channels</p> <div data-bbox="722 321 1820 459"> <div data-bbox="722 321 1073 459"> <p>Finely scaled frequency resolution in 20 channels gives the Fine-scale Noise Canceler (FNC) the most precise background noise reduction capability. Noise reduction depends on the relative amounts of speech and noise in each channel (graph 1) and the importance of each</p> </div> <div data-bbox="1098 321 1449 459"> <p>channel for speech understanding. The FNC in Perseo ensures maximum comfort and ease of listening. Annoying background sounds are reduced for comfort while the important sounds of speech are preserved for stress-free listening. The maximum amount of noise</p> </div> <div data-bbox="1474 321 1820 459"> <p>reduction can be fitter-programmed to match the individual wearer's tastes and requirements. Thanks to the dual-base automatic, separate settings can be made for relatively quiet and noisy situations, ensuring optimal comfort in multiple situations (graph 2).</p> </div> </div>  <p data-bbox="583 841 917 873">Anderson '783 at 2:20–39</p> <p data-bbox="632 878 1988 1276">It is an object of the present invention to provide a new and useful auditory aid for hearing-impaired persons (i.e., those having certain residual hearing) by removing audio signal enhancement functions from the earpiece and placing them in an RPU. Use of an RPU provides several advantages over systems that attempt to place all system capabilities within the earpiece. The RPU approach allows a simple earpiece design comprising a miniature low-power wireless transceiver, microphone and speaker. Note that the speaker is also known in the hearing aid art as a "receiver," but the term "speaker" is used here to avoid confusion (similarly, the term "talker," not "speaker," is used to describe a human producing vocal sounds). The resulting earpiece is extremely small, can be hidden from view in the ear canal if desired, and allows complete freedom of movement when a primary two-way wireless link to the RPU is used. Processing for all major system capabilities, such as amplification and other forms of signal enhancement, takes place in the RPU where size and power constraints are relaxed, leading to a cost-effective design.</p> <p data-bbox="583 1333 974 1365">Anderson '783 at 22:61–23:35</p> <p data-bbox="632 1370 1988 1476">Details of normal system operation: During normal operation of the hearing aid with wireless remote processor, sounds from the ambient environment are picked up by a microphone 12 in the earpiece 10 and transmitted with other information (e.g., data representing the earpiece battery 860 voltage as obtained by</p>

Claims	Prior Art
	<p>an A/D 865 converter and transmitted under control of the earpiece ASIC 845 in the auxiliary bit of many successive replies via the reply data register 825) over a primary two-way wireless link 17 to the RPU 16, where the audio signals are enhanced according to the user's needs by the RPU DSP 948. Signal enhancement is achieved via any of a number of techniques well known in the art that use the capabilities of the general-purpose RPU DSP 948. The enhanced audio signals may be combined with other information (e.g., a synthesized voice, generated within the DSP 948, warning that the aforementioned earpiece battery voltage is low) and transmitted from the RPU 16 over the primary wireless link 17 to the earpiece 10, where they are converted by a speaker 15 to sounds that can only be heard by the user 11. Art auxiliary speaker can be supplied in the RPU 16 as a peripheral 950 device that allows people other than the user to hear the same sounds as those produced by the earpiece speaker 15, if desired (as, e.g., when a person other than the user wishes to use the system's optional cellular telephone capabilities). A different, optional secondary two-way wireless link (shown as part of the RPU in FIG. 9 as secondary wireless link circuitry 944, or connected to the RPU by a wired or other, e.g. infrared, wireless link 18 and shown as optional secondary wireless link circuitry 19 in FIG. 1) can be used for communication between the RPU and a cellular telephone system or other sources of information. An RPU keyboard suitable for covert data entry 72, or voice recognition capabilities in the RPU (implemented, e.g., using any of a number of techniques well known in the art that use the general-purpose RPU DSP 948), can be used to control hearing aid parameters (such as amplification level set by the RPU DSP 948) and telephone dialing functions (e.g., automatic dialing of a predefined telephone number using information stored in the RPU DSP 948) in a manner that can be made imperceptible to a casual observer if desired.</p> <p>Anderson '783 at 27:4–24</p> <p>Hearing test: The present invention allows hearing tests to be performed via the earpiece normally worn by the user, without the need for additional expensive test equipment. All signal generation capabilities necessary to perform hearing tests are available in the hearing aid system comprising an earpiece 10 and RPU 16. Specifically, the RPU DSP 948 can be used to generate tones of varying frequencies and amplitudes, as well as other signals known to be useful for hearing test purposes (e.g., nonsense syllables) that are converted to acoustic waves by the earpiece. Such tones are used to perform a hearing test that determines appropriate gain vs. frequency parameters for a program stored in the RPU DSP 948 that performs signal enhancement to compensate for the user's hearing loss. Note that the hearing test program that controls the RPU DSP 948 during the hearing test can be temporarily stored in the RPU DSP 948 for the duration of the test, then deleted upon completion of the test to allow re-use of RPU DSP 948 memory resources during normal operation. The hearing test program may be loaded into the RPU DSP 948</p>

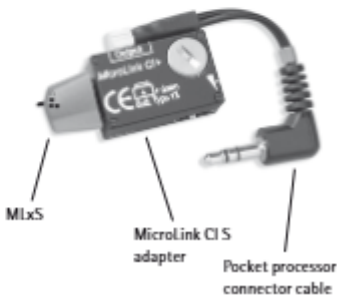
Claims	Prior Art
	<p>through the secondary wireless link 944 or a wired peripheral link 950.</p> <p>A POSITA would have been motivated to combine the teachings of the DHA system with Anderson to meet the claim limitations with a reasonable expectation of success for the reasons described above with respect to claim limitations [1b] and [1c].</p> <p>To the extent that the DHA system and Anderson not disclose any elements of this claim limitation, it would have been obvious to a person having ordinary skill in the art to combine, with the teachings of the DHA system, those elements, as known to them or as disclosed by other prior art identified in the cover pleading.</p>
<p>[1e] the mobile communication device configured to transmit the customized audio output signals to the pair of ear buds,</p>	<p>DHA systems, alone or in view of Anderson, expressly or inherently discloses “the mobile communication device configured to transmit the customized audio output signals to the pair of ear buds.”</p> <p>Phonak - SmartLink SX - FM Systems Manual at p. 3</p> <p>SmartLink SX is the first communication device to incorporate</p> <ul style="list-style-type: none"> • Advanced FM features • Hearing Instrument remote control • A Bluetooth link for cell phone use <p>Phonak - SmartLink SX - FM Systems Manual at p. 5</p> <p>Smart FM transmission</p> <p>The proven Phonak FM technology clearly transmits all incoming auditory signals from SmartLink to the receiver attached to the user's hearing instrument.</p> <p>Phonak - SmartLink SX User Guide at p. 6</p> <p>This unique combination of an FM transmitter, mobile phone link and hearing instrument remote control is the gateway to a completely new hearing experience.</p>

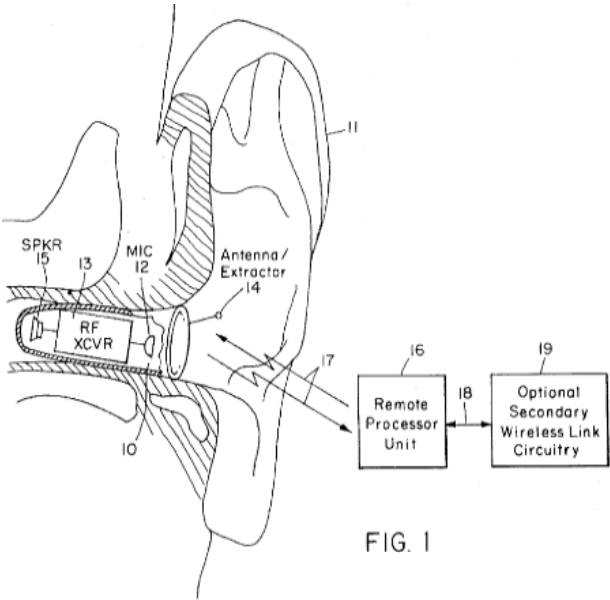
Claims	Prior Art
	<p>Phonak - SmartLink SX User Guide at p. 32</p>  <p>Phonak - SmartLink SX User Guide at p. 39</p> <p><u>Answering a call</u></p> <ul style="list-style-type: none"> • Speak directly into the SmartLink. You hear the caller's voice in both hearing instruments. <p>Phonak - SmartLink SX User Guide at p. 58</p> <p>The SmartLink is used by hearing-impaired people. To operate the SmartLink, BTE hearing instruments with MicroLink receivers are required.</p> <p>Phonak – Microlink User Guide at Page 2</p> <p>The desired signal – for example, a speaker's voice or the sound from a TV – is picked up by the microphone which is placed near to the source of the sound. The transmitter sends an FM signal to the receiver, which immediately relays this to the hearing instrument for you to hear – crystal clear!</p>

Claims	Prior Art
	<p>Phonak – Online Interview re Perseo (April 7, 2003) at p. .5</p> <p><u>AO/Beck</u>: Phonak has been a leader in FM technology-is this available for the Perseo product line?</p> <p><u>Dave Fabry/Phonak</u>: Yes, FM can be an integrated solution for the Perseo BTEs and is compatible with all Phonak transmitters.</p> <p>Phonak – Claro 211 Product Information at Page 2</p> <p><u>Principle of FM operation</u></p> <p>The transmitter unit is worn by the person speaking, placed in front of the sound source or even connected directly to TV or HiFi. The transmitter unit collects the desired signal and transmits it via radio waves directly to the FM module integrated in the hearing instrument. The incoming FM signal received by the hearing instrument is processed either in accordance with the audiological setting of the instrument, providing consistent sound quality or according to an individual FM program.</p> <p>Anderson '783 at FIG. 1</p>  <p>Anderson '783 at 22:61–23:35</p> <p>Details of normal system operation: During normal operation of the hearing aid with wireless remote processor, sounds from the ambient environment are picked up by a microphone 12 in the earpiece 10 and</p>

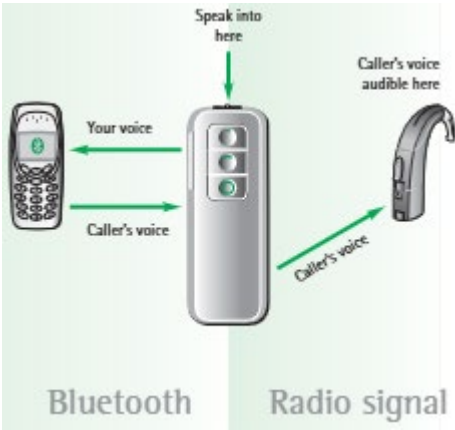
Claims	Prior Art
	<p>transmitted with other information (e.g., data representing the earpiece battery 860 voltage as obtained by an A/D 865 converter and transmitted under control of the earpiece ASIC 845 in the auxiliary bit of many successive replies via the reply data register 825) over a primary two-way wireless link 17 to the RPU 16, where the audio signals are enhanced according to the user's needs by the RPU DSP 948. Signal enhancement is achieved via any of a number of techniques well known in the art that use the capabilities of the general-purpose RPU DSP 948. The enhanced audio signals may be combined with other information (e.g., a synthesized voice, generated within the DSP 948, warning that the aforementioned earpiece battery voltage is low) and transmitted from the RPU 16 over the primary wireless link 17 to the earpiece 10, where they are converted by a speaker 15 to sounds that can only be heard by the user 11. Art auxiliary speaker can be supplied in the RPU 16 as a peripheral 950 device that allows people other than the user to hear the same sounds as those produced by the earpiece speaker 15, if desired (as, e.g., when a person other than the user wishes to use the system's optional cellular telephone capabilities). A different, optional secondary two-way wireless link (shown as part of the RPU in FIG. 9 as secondary wireless link circuitry 944, or connected to the RPU by a wired or other, e.g. infrared, wireless link 18 and shown as optional secondary wireless link circuitry 19 in FIG. 1) can be used for communication between the RPU and a cellular telephone system or other sources of information. An RPU keyboard suitable for covert data entry 72, or voice recognition capabilities in the RPU (implemented, e.g., using any of a number of techniques well known in the art that use the general-purpose RPU DSP 948), can be used to control hearing aid parameters (such as amplification level set by the RPU DSP 948) and telephone dialing functions (e.g., automatic dialing of a predefined telephone number using information stored in the RPU DSP 948) in a manner that can be made imperceptible to a casual observer if desired.</p> <p>A POSITA would have been motivated to combine the teachings of the DHA system with Anderson to meet the claim limitations with a reasonable expectation of success for the reasons described above with respect to claim limitation [1b] and [1c].</p> <p>To the extent that the DHA system and Anderson not disclose any elements of this claim limitation, it would have been obvious to a person having ordinary skill in the art to combine, with the teachings of the DHA system, those elements, as known to them or as disclosed by other prior art identified in the cover pleading.</p>
[1f] wherein each of the ear buds is configured to receive at least a portion of the customized audio output signals transmitted by the mobile communication device, and	DHA systems, alone or in view of Anderson, expressly or inherently discloses “wherein each of the ear buds is configured to receive at least a portion of the customized audio output signals transmitted by the mobile communication device.”

Claims	Prior Art
	<p>Phonak - SmartLink SX - FM Systems Manual at p. 5</p> <p>Smart FM transmission The proven Phonak FM technology clearly transmits all incoming auditory signals from SmartLink to the receiver attached to the user's hearing instrument.</p> <p>Phonak - SmartLink SX User Guide at p. 6</p> <p>This unique combination of an FM transmitter, mobile phone link and hearing instrument remote control is the gateway to a completely new hearing experience.</p> <p>Phonak - SmartLink SX User Guide at p. 32</p>  <p>Phonak – Microlink User Guide at Page 2</p> <p>The desired signal – for example, a speaker's voice or the sound from a TV – is picked up by the microphone which is placed near to the source of the sound. The transmitter sends an FM signal to the receiver, which immediately relays this to the hearing instrument for you to hear – crystal clear!</p>

Claims	Prior Art
	<p>Phonak – Claro 211 Product Information at Page 2</p> <p>Principle of FM operation</p> <p>The transmitter unit is worn by the person speaking, placed in front of the sound source or even connected directly to TV or HiFi. The transmitter unit collects the desired signal and transmits it via radio waves directly to the FM module integrated in the hearing instrument. The incoming FM signal received by the hearing instrument is processed either in accordance with the audiological setting of the instrument, providing consistent sound quality or according to an individual FM program.</p> <p>Phonak MicroLink MLxS Product Information at Page 1</p> <p>The MicroLink MLxS is the first miniaturized FM receiver with a digital frequency synthesizer for full multi-frequency flexibility. The receiver features fully automatic, non-contact frequency synchronization with either a WallPilot or the Campus S.</p> <p>The MicroLink MLxS is compatible with all Phonak FM transmitters and the WallPilot, and is recommended for use with this equipment.</p> <div style="display: flex; justify-content: space-around; align-items: flex-start;"> <div style="text-align: center;"> <p>MLxS with BTE instrument</p>  </div> <div style="text-align: center;"> <p>MLxS with adapter for CI pocket processor</p>  </div> </div>

Claims	Prior Art
	<p data-bbox="583 175 894 207">Anderson '783 at FIG. 1</p>  <p data-bbox="995 748 1062 776">FIG. 1</p> <p data-bbox="583 846 919 878">Anderson '783 at 3:13–18</p> <p data-bbox="632 883 1969 1024">FIG. 1 is a block diagram system overview illustrating how an earpiece (comprising a microphone, RF transceiver, dual-use antenna/extractor and speaker) is worn in the ear, communicates via a primary two-way RF link with an RPU, and also communicates via an optional secondary wireless link to a telephone system.</p> <p data-bbox="583 1081 945 1114">Anderson '783 at 3:52–4:14</p> <p data-bbox="632 1118 1986 1442">One preferred embodiment of the invention herein described is shown in block diagram form in FIG. 1. In this figure, the earpiece 10 shown worn in the ear 11 uses a standard completely-in-the-canal (CIC) housing, well known in the art, although many other earpiece housing types known in the art (e.g., behind-the-ear, or BTE) can be used. The earpiece shown in FIG. 1 comprises a microphone 12, RF transceiver 13 with an antenna that doubles as an earpiece extractor 14, and speaker 15. Although the dual-use antenna/extractor arrangement shown in FIG. 1 is often desirable, a separate antenna or antennae can be used in conjunction with many different types of earpiece extractors as appropriate for a particular application. Antennae and extractors can be disguised as ear hair or jewelry (e.g., earrings), or conductive filaments can be permanently implanted in the ear cartilage to act as antennae using a process similar to</p>

Claims	Prior Art
	<p>ear piercing. The earpiece 10 communicates with a remote processor unit (RPU) 16 via a primary two-way RF link 17, although many other wireless link media (e.g., ultrasonic or infrared) may be used instead. The RPU is typically worn under clothing (e.g., carried in a pocket or purse), but may also be worn in plain sight (e.g., on a belt) if desired. The RPU 16 may be connected (via wired or wireless means 18) to optional secondary wireless link circuitry 19 that allows wireless communication between the RPU and other sources of information (e.g., the general subscriber telephone network) via a secondary wireless link. Note that the optional secondary wireless link circuitry 19 may or may not be contained within the RPU case 70 of FIG. 7.</p> <p>Anderson '783 at 4:30–33</p> <p>Although the system described here is monaural, it will be clear to those skilled in the art that a binaural (stereo) system can be constructed based on the same principles as the monaural system.</p> <p>Anderson '783 at 22:61–23:35</p> <p>Details of normal system operation: During normal operation of the hearing aid with wireless remote processor, sounds from the ambient environment are picked up by a microphone 12 in the earpiece 10 and transmitted with other information (e.g., data representing the earpiece battery 860 voltage as obtained by an A/D 865 converter and transmitted under control of the earpiece ASIC 845 in the auxiliary bit of many successive replies via the reply data register 825) over a primary two-way wireless link 17 to the RPU 16, where the audio signals are enhanced according to the user's needs by the RPU DSP 948. Signal enhancement is achieved via any of a number of techniques well known in the art that use the capabilities of the general-purpose RPU DSP 948. The enhanced audio signals may be combined with other information (e.g., a synthesized voice, generated within the DSP 948, warning that the aforementioned earpiece battery voltage is low) and transmitted from the RPU 16 over the primary wireless link 17 to the earpiece 10, where they are converted by a speaker 15 to sounds that can only be heard by the user 11. Art auxiliary speaker can be supplied in the RPU 16 as a peripheral 950 device that allows people other than the user to hear the same sounds as those produced by the earpiece speaker 15, if desired (as, e.g., when a person other than the user wishes to use the system's optional cellular telephone capabilities). A different, optional secondary two-way wireless link (shown as part of the RPU in FIG. 9 as secondary wireless link circuitry 944, or connected to the RPU by a wired or other, e.g. infrared, wireless link 18 and shown as optional secondary wireless link circuitry 19 in FIG. 1) can be used for communication between the RPU and a cellular telephone system or other sources of information. An RPU keyboard suitable for covert data entry 72, or voice recognition capabilities in the RPU (implemented, e.g., using any of a number of techniques well known in the art that use the general-purpose RPU DSP 948), can be used to control</p>

Claims	Prior Art
	<p>hearing aid parameters (such as amplification level set by the RPU DSP 948) and telephone dialing functions (e.g., automatic dialing of a predefined telephone number using information stored in the RPU DSP 948) in a manner that can be made imperceptible to a casual observer if desired.</p> <p>A POSITA would have been motivated to combine the teachings of the DHA system with Anderson to meet the claim limitations with a reasonable expectation of success for the reasons described above with respect to claim limitation [1b] and [1c].</p> <p>To the extent that the DHA system, alone or in view of Anderson, does not disclose any elements of this claim limitation, it would have been obvious to a person having ordinary skill in the art to combine, with the teachings of DHA systems, those elements, as known to them or as disclosed by other prior art identified in the cover pleading.</p>
<p>[1g] wherein each of the ear buds is configured to generate at least a portion of the customized audio output for the user in accordance with at least the portion of the customized audio output signals being received.</p>	<p>DHA systems, alone or in view of Anderson, expressly or inherently discloses “wherein each of the ear buds is configured to generate at least a portion of the customized audio output for the user in accordance with at least the portion of the customized audio output signals being received.”</p> <p>Phonak - SmartLink SX User Guide at p. 32</p>  <p>Phonak - SmartLink SX User Guide at p. 39</p> <p><u>Answering a call</u></p> <ul style="list-style-type: none"> • Speak directly into the SmartLink. You hear the caller's voice in both hearing instruments.

Claims	Prior Art
	<p data-bbox="583 175 1150 207">Phonak - SmartLink SX User Guide at p. 58</p> <p data-bbox="638 228 1062 298">The user benefit of FM is an improved signal-to-noise ratio, which leads to an optimized speech understanding even in very difficult listening situations.</p> <p data-bbox="583 326 1062 358">Phonak – Perseo Brochure at Page 13</p> <div data-bbox="627 415 1482 456"> <h2>Enhanced sound quality and Speech-to-Noise Ratio</h2>  </div> <p data-bbox="726 509 1062 532">Clarity Component Enhancement – DPP²</p> <p data-bbox="726 550 1719 764">Speech clarity depends on both its acoustic properties and the listener's auditory perception. Conventional spectral enhancement systems consider only the acoustics of sound, not how it will be perceived by the listener. Perseo goes well beyond this with DPP², the new generation of Digital Perception Processing. DPP² with Clarity Component Enhancement is individualized for each wearer. It analyzes the acoustic and the psychoacoustic characteristics of incoming sound in 20 critical bands to identify accurately the components that contribute to clarity. It then selectively amplifies them to optimize clarity for the wearer.</p> <p data-bbox="583 786 1117 818">Phonak – Microlink User Guide at Page 2</p> <p data-bbox="638 839 1157 993">The desired signal – for example, a speaker's voice or the sound from a TV – is picked up by the microphone which is placed near to the source of the sound. The transmitter sends an FM signal to the receiver, which immediately relays this to the hearing instrument for you to hear – crystal clear!</p> <p data-bbox="583 1024 1230 1057">Phonak – Claro 211 Product Information at Page 2</p> <p data-bbox="638 1083 877 1105">Principle of FM operation</p> <p data-bbox="638 1131 1167 1274">The transmitter unit is worn by the person speaking, placed in front of the sound source or even connected directly to TV or HiFi. The transmitter unit collects the desired signal and transmits it via radio waves directly to the FM module integrated in the hearing instrument. The incoming FM signal received by the hearing instrument is processed either in accordance with the audiological setting of the instrument, providing consistent sound quality or according to an individual FM program.</p>

Claims

Prior Art

Phonak – Collections 2002 Product List at Page 7



Adaptive digital AudioZoom (dAZ)



Fine-scale Noise Canceler (FNC)

Active and adaptive separation of speech from noise
The Speech-to-Noise Ratio (SNR) can only be optimised if the noise is spatially separated from the speech signal with directional microphone technology. As a pioneer in multi-microphone technologies, Phonak introduced AudioZoom, a feature which significantly improves listening and understanding in noisy situations, enhancing user satisfaction. Now, Phonak goes a step further in the refinement of multi-microphone technology with Adaptive digital AudioZoom: it seeks, finds, and suppresses noise automatically.

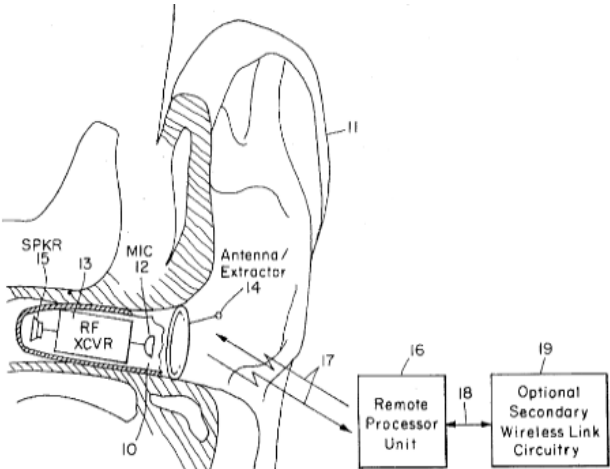
Fully automatic, fully adaptive
With its unique spatial processing system, Phonak's Adaptive digital AudioZoom can detect noise not only from directly behind, but anywhere in the near plane. Whether the hearing-impaired person is moving or the location of the noise source changes, the directional characteristic is constantly adapted to the listening situation. In this manner, the Adaptive digital AudioZoom automatically establishes the best possible SNR.

Twenty-band noise control
Claro constantly monitors the Speech-to-Noise Ratio (SNR) in all 20 bands. If a poor SNR is measured in certain bands, the Fine-scale Noise Canceler will reduce gain only in those bands. The result of this precision noise control is a selective and intelligent improvement of the SNR with only a marginal influence on the speech signal.

Intelligent
The magnitude of gain reduction is not constant. In frequency bands with higher SNRs and in bands which are important for speech intelligibility, gain is reduced only slightly. This intelligent control algorithm reduces the masking and distractive effects of background noise. Important speech information is preserved.

Precise and selective
All frequency bands are separately analysed and controlled. Gain is reduced selectively only in those bands where noise is present. This high-resolution analysis and gain control makes it possible to suppress noise with extremely high accuracy. Neighboring frequency bands without noise remain unaffected.

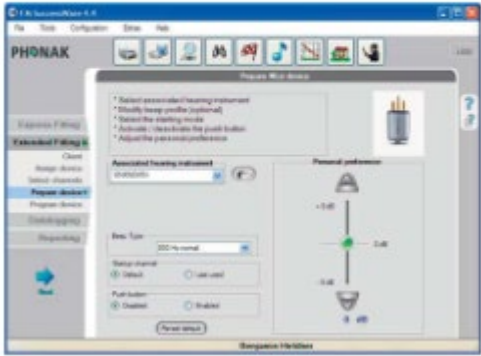
Claims	Prior Art
	<p data-bbox="583 175 1608 212">Lewis – Speech perception in noise: directional microphones (2004) at Page 431</p> <p data-bbox="636 228 1016 818"> All subjects were fit with digital Phonak Claro 311 dAZ BTE hearing aids bilaterally. All earmolds had select-a-vent (SAV) venting and #13 or 3 mm horn tubing. In addition to the hearing aids, subjects were fit with Phonak Microlink ML8 FM receivers bilaterally. These FM receivers attach to the bottom of a BTE hearing aid and may be utilized in either the “FM only” mode, which attenuates the hearing aid microphone by 20 dB, or in the FM plus hearing aid mode, which allows for FM input and input of environmental sounds via the hearing aid microphones simultaneously at the same output level. The Phonak TX3 HandyMic FM transmitter served as the FM transmitter. This transmitter has three microphone options: (1) Wide Angle, which picks up sounds arriving from all directions around the transmitter microphone equally; (2) Zoom, which provides reduced amplification to signals arriving from the rear (cardioid); and (3) SuperZoom, which provides reduced amplification to signals originating from the rear and the sides (hypercardioid). The </p> <p data-bbox="583 846 1115 878">Oticon – DigiFocus Description at Page 1</p> <p data-bbox="636 899 1314 1062"> Inside each DigiFocus, the digital chips have the same processing power as a fast laptop computer. This power is used to analyse the incoming sound signal, and speech is extracted through a unique system called ASA - Adaptive Speech Alignment. ASA always makes speech audible, amplifying it as much as necessary without making it uncomfortably loud or distorted. Sounds are reproduced very naturally - making speech far easier to understand than with any other hearing instrument on the market. </p> <p data-bbox="583 1089 1220 1122">Oticon – DigiFocus White Paper (1996) at Page 3</p> <p data-bbox="636 1138 1299 1349"> After the A/D-conversion, the resulting digital code can be mathematically manipulated, or processed (as it is normally expressed when talking about sound). This is where the real benefits of digital technology come in. The amount and types of processing that are possible with digital technology far surpass the capabilities of analogue processing. Once the processing is complete, the digital signal is converted back to analogue form by the D/A-converter. It is now ready for transmission to the small receiver inside the hearing instrument, which makes the sound audible. □ </p>

Claims	Prior Art
	<p data-bbox="583 175 894 207">Anderson '783 at FIG. 1</p>  <p data-bbox="995 748 1062 781">FIG. 1</p> <p data-bbox="583 846 945 878">Anderson '783 at 3:52–4:14</p> <p data-bbox="627 883 1992 1463">One preferred embodiment of the invention herein described is shown in block diagram form in FIG. 1. In this figure, the earpiece 10 shown worn in the ear 11 uses a standard completely-in-the-canal (CIC) housing, well known in the art, although many other earpiece housing types known in the art (e.g., behind-the-ear, or BTE) can be used. The earpiece shown in FIG. 1 comprises a microphone 12, RF transceiver 13 with an antenna that doubles as an earpiece extractor 14, and speaker 15. Although the dual-use antenna/extractor arrangement shown in FIG. 1 is often desirable, a separate antenna or antennae can be used in conjunction with many different types of earpiece extractors as appropriate for a particular application. Antennae and extractors can be disguised as ear hair or jewelry (e.g., earrings), or conductive filaments can be permanently implanted in the ear cartilage to act as antennae using a process similar to ear piercing. The earpiece 10 communicates with a remote processor unit (RPU) 16 via a primary two-way RF link 17, although many other wireless link media (e.g., ultrasonic or infrared) may be used instead. The RPU is typically worn under clothing (e.g., carried in a pocket or purse), but may also be worn in plain sight (e.g., on a belt) if desired. The RPU 16 may be connected (via wired or wireless means 18) to optional secondary wireless link circuitry 19 that allows wireless communication between the RPU and other sources of information (e.g., the general subscriber telephone network) via a secondary wireless link. Note that the optional secondary wireless link circuitry 19 may or may not be contained within the</p>

Claims	Prior Art
	<p>RPU case 70 of FIG. 7.</p> <p>Anderson '783 at 4:30–33 Although the system described here is monaural, it will be clear to those skilled in the art that a binaural (stereo) system can be constructed based on the same principles as the monaural system.</p> <p>Anderson '783 at 22:11–61 Stand-alone earpiece operation: If an earpiece cannot communicate with an RPU for any reason, stand-alone earpiece operation as a CIC hearing aid is generally desirable. Examples include occasions where the RPU must be turned off by pushing the RPU power-off key 78 (e.g., situations where operation of wireless transmitters is prohibited such as aircraft takeoffs and classified military briefings), RPU battery or other RPU electronics failure, misplacement of an RPU, and primary wireless link 17 interference. Note that a microphone located in the earpiece 10 is preferred to a body-worn microphone, since only an earpiece microphone 12 can be used to provide the desired default CIC hearing aid function (earpiece microphones also provide superior audio localization capabilities as the microphone moves with the user's head, and provide reduced wind and clothing noise compared with body-worn microphones). Parameters that control the stand-alone earpiece operation are stored in the earpiece CIC default parameter storage 870, which is accessed by the ASIC 845 to implement the desired stand-alone functions as described earlier. CIC default parameters include, but are not limited to, settings for the nominal and maximum allowable earpiece volume. Other information contained in the default parameter storage 870 includes, but is not limited to, the communication address and L/R (left/right) selection bit. To ensure that the CIC default parameter storage information cannot be changed inadvertently, changes are only allowed when the earpiece is not in use as, for example, when the earpiece battery 860 is being charged. A battery charging detector 840, that is sensitive to the high-level magnetic induction fields commonly used for earpiece battery charging, indicates to the ASIC 845 that the default parameter storage 870 values may be changed upon receipt of an appropriate command from the RPU. Specifically, to change information in the earpiece CIC default parameter storage 870, a user first activates the earpiece battery charging system, e.g., by placing the earpiece in a recharger/storage tray well known in the art. The RPU keyboard 946 is then used to enter a command to the DSP 948 that indicates the user's desire to change information in the default parameter storage 870, followed by the new information to be stored. The DSP 948 causes the RPU to transmit an interrogation with a reserved address code to which all earpieces being charged within the communication range of the RPU are receptive. The information that is to ultimately be loaded into the earpiece default parameter storage 870 is then transmitted by the RPU using the auxiliary bit field of many successive interrogations, said interrogations are received by the earpiece, and the desired</p>

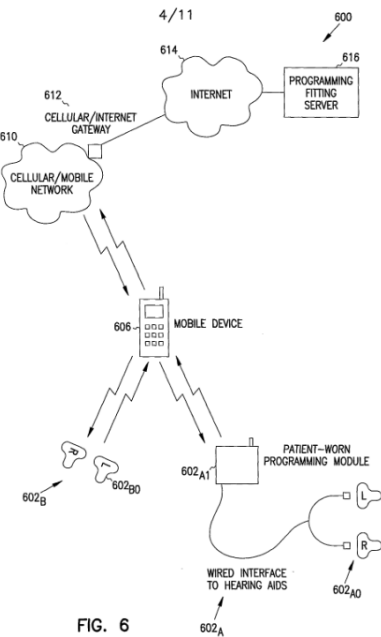
Claims	Prior Art
	<p>information is loaded into the CIC default parameter storage 870 under control of the ASIC 845.</p> <p>Anderson '783 at 22:61–23:35</p> <p>Details of normal system operation: During normal operation of the hearing aid with wireless remote processor, sounds from the ambient environment are picked up by a microphone 12 in the earpiece 10 and transmitted with other information (e.g., data representing the earpiece battery 860 voltage as obtained by an A/D 865 converter and transmitted under control of the earpiece ASIC 845 in the auxiliary bit of many successive replies via the reply data register 825) over a primary two-way wireless link 17 to the RPU 16, where the audio signals are enhanced according to the user's needs by the RPU DSP 948. Signal enhancement is achieved via any of a number of techniques well known in the art that use the capabilities of the general-purpose RPU DSP 948. The enhanced audio signals may be combined with other information (e.g., a synthesized voice, generated within the DSP 948, warning that the aforementioned earpiece battery voltage is low) and transmitted from the RPU 16 over the primary wireless link 17 to the earpiece 10, where they are converted by a speaker 15 to sounds that can only be heard by the user 11. Art auxiliary speaker can be supplied in the RPU 16 as a peripheral 950 device that allows people other than the user to hear the same sounds as those produced by the earpiece speaker 15, if desired (as, e.g., when a person other than the user wishes to use the system's optional cellular telephone capabilities). A different, optional secondary two-way wireless link (shown as part of the RPU in FIG. 9 as secondary wireless link circuitry 944, or connected to the RPU by a wired or other, e.g. infrared, wireless link 18 and shown as optional secondary wireless link circuitry 19 in FIG. 1) can be used for communication between the RPU and a cellular telephone system or other sources of information. An RPU keyboard suitable for covert data entry 72, or voice recognition capabilities in the RPU (implemented, e.g., using any of a number of techniques well known in the art that use the general-purpose RPU DSP 948), can be used to control hearing aid parameters (such as amplification level set by the RPU DSP 948) and telephone dialing functions (e.g., automatic dialing of a predefined telephone number using information stored in the RPU DSP 948) in a manner that can be made imperceptible to a casual observer if desired.</p> <p>A POSITA would have been motivated to combine the teachings of the DHA system with Anderson to meet the claim limitations with a reasonable expectation of success for the reasons described above with respect to claim limitation [1b] and [1c].</p> <p>To the extent that the DHA system, alone or in view of Anderson, does not disclose any elements of this claim limitation, it would have been obvious to a person having ordinary skill in the art to combine, with the teachings of DHA systems, those elements, as known to them or as disclosed by other prior art identified in the cover pleading.</p>

Claims	Prior Art
<p>2. A personal audio system as recited in claim 1, wherein the mobile communication device is configured to establish the hearing profile for the user.</p>	<p>DHA systems, alone or in view of Anderson, expressly or inherently discloses “[a] personal audio system as recited in claim 1, wherein the mobile communication device is configured to establish the hearing profile for the user.”</p> <p>Phonak – Collections 2002 Product List at Page 33</p> <div data-bbox="630 354 1302 950">   <p>PFG – new perspectives in fitting success</p> <p>QuickGuide The versatile, easy-to-use online tour through the PFG fitting process with access to all optional fitting tools.</p> <p>Instrument selection wizard This online assistant recommends the product that best matches the customer's needs.</p> <p>Initial SoundCheck This is an ideal tool for fast and simple fine tuning of sound quality. Three controls provide convenient gain and balance settings in the high- and low-frequency ranges.</p> <p>Automatic fine tuning Based on the customer's hearing perception, the software automatically computes the most suitable modifications.</p> <p>Manual fine tuning A high degree of fitting precision can be attained with this broad spectrum of parameters in an optimised display format.</p> <p>Phonak Sound The degree of fitting success can be dramatically improved by using these real-life sound clips. All sound examples can be saved on a computer as MP3 files; they are also available on a sound CD.</p> <p>Proven PFG benefits</p> <ul style="list-style-type: none"> • One program for all Phonak hearing computers • NOAH compatible or standalone version • Detailed help system with in-depth scientific and technical information • PFG training demo </div> <p>Anderson '783 at 27:4–24</p> <p>Hearing test: The present invention allows hearing tests to be performed via the earpiece normally worn by the user, without the need for additional expensive test equipment. All signal generation capabilities necessary to perform hearing tests are available in the hearing aid system comprising an earpiece 10 and RPU 16. Specifically, the RPU DSP 948 can be used to generate tones of varying frequencies and amplitudes, as well as other signals known to be useful for hearing test purposes (e.g., nonsense syllables) that are converted to acoustic waves by the earpiece. Such tones are used to perform a hearing test that determines appropriate gain vs. frequency parameters for a program stored in the RPU DSP 948 that performs signal enhancement to compensate for the user's hearing loss. Note that the hearing test program that controls the RPU DSP 948 during the hearing test can be temporarily stored in the RPU DSP 948 for the duration of the test, then deleted upon completion of the test to allow re-use of RPU DSP 948 memory resources during normal operation. The hearing test program may be loaded into the RPU DSP 948</p>

Claims	Prior Art
	<p>through the secondary wireless link 944 or a wired peripheral link 950.</p> <p>A POSITA would have been motivated to combine the teachings of the DHA system with Anderson to meet the claim limitations with a reasonable expectation of success for the reasons described above with respect to claim limitation [1b].</p> <p>To the extent that the DHA system, alone or in view of Anderson, does not disclose any elements of this claim limitation, it would have been obvious to a person having ordinary skill in the art to combine, with the teachings of DHA systems, those elements, as known to them or as disclosed by other prior art identified in the cover pleading.</p>
<p>3. A personal audio system as recited in claim 2, wherein the hearing profile for the user is established by the mobile communication device at least via generating audio tones and outputting such to one or more of the ear buds and receiving a user indication corresponding to the user's ability to hear the generated audio tones.</p>	<p>DHA systems, alone or in view of Anderson or Eaton, expressly or inherently discloses “[a] personal audio system as recited in claim 2, wherein the hearing profile for the user is established by the mobile communication device at least via generating audio tones and outputting such to one or more of the ear buds and receiving a user indication corresponding to the user's ability to hear the generated audio tones.”</p> <p>Phonak – Fitting Guide for MLxi and MLxi Baha at Page 4</p> <p>4.2 Advanced fitting of MLxi using FM SuccessWare</p> <p>The extended mode of FM SuccessWare enables you to fine-tune some parameters of the MLxi according to the needs of your customers.</p> 

Claims	Prior Art
	<p>Phonak – Collections 2002 Product List at Page 33</p> <div data-bbox="636 256 802 337"> <p>PFG – new perspectives in fitting success</p> </div> <div data-bbox="636 386 884 570">  </div> <div data-bbox="636 591 884 773">  </div> <div data-bbox="945 227 1302 820"> <p>QuickGuide The versatile, easy-to-use online tour through the PFG fitting process with access to all optional fitting tools.</p> <p>Instrument selection wizard This online assistant recommends the product that best matches the customer's needs.</p> <p>Initial SoundCheck This is an ideal tool for fast and simple fine tuning of sound quality. Three controls provide convenient gain and balance settings in the high- and low-frequency ranges.</p> <p>Automatic fine tuning Based on the customer's hearing perception, the software automatically computes the most suitable modifications.</p> <p>Manual fine tuning A high degree of fitting precision can be attained with this broad spectrum of parameters in an optimised display format.</p> <p>Phonak Sound The degree of fitting success can be dramatically improved by using these real-life sound clips. All sound examples can be saved on a computer as MP3 files; they are also available on a sound CD.</p> <p>Proven PFG benefits</p> <ul style="list-style-type: none"> • One program for all Phonak hearing computers • NOAH compatible or standalone version • Detailed help system with in-depth scientific and technical information • PFG training demo </div> <p>Anderson '783 at 27:4–24</p> <p>Hearing test: The present invention allows hearing tests to be performed via the earpiece normally worn by the user, without the need for additional expensive test equipment. All signal generation capabilities necessary to perform hearing tests are available in the hearing aid system comprising an earpiece 10 and RPU 16. Specifically, the RPU DSP 948 can be used to generate tones of varying frequencies and amplitudes, as well as other signals known to be useful for hearing test purposes (e.g., nonsense syllables) that are converted to acoustic waves by the earpiece. Such tones are used to perform a hearing test that determines appropriate gain vs. frequency parameters for a program stored in the RPU DSP 948 that performs signal enhancement to compensate for the user's hearing loss. Note that the hearing test program that controls the RPU DSP 948 during the hearing test can be temporarily stored in the RPU DSP 948 for the duration of the test, then deleted upon completion of the test to allow re-use of RPU DSP 948 memory resources during normal operation. The hearing test program may be loaded into the RPU DSP 948 through the secondary wireless link 944 or a wired peripheral link 950.</p>

Claims	Prior Art
	<p data-bbox="583 175 1980 284">To the extent that Anderson '783 does not disclose any elements of this claim limitation, it would have been known to a person having ordinary skill in the art at the time of the invention that executing software to test a user's hearing is a common part of fitting a user for a hearing aid. For example, Eaton '458 teaches:</p> <p data-bbox="583 300 869 332">Eaton '458 at 11:6–15</p> <p data-bbox="630 349 1980 641">Figure 6 is a pictorial diagram illustrating a system in accordance with one embodiment. Figure 6 contains elements similar to those discussed in Figure 1. For these elements, the aforementioned discussion related to them is incorporated here in full. The system 600 includes hearing aid systems 602_A and 602_B. The hearing aid system 602_A comprises a programming module 602_{A1} and a hearing aid device 602_{A0}. The hearing aid system 602_B comprises the hearing aid device 602_{B0}. Whereas a mobile device 606 communicates with the hearing aid device 602_{A0} of the hearing aid system 602_A through the programming module 602_{A1}, the mobile device 606 communicates directly with the hearing aid device 602_{B0} of the hearing aid system 602_B.</p> <p data-bbox="583 657 884 690">Eaton '458 at 11:16–28</p> <p data-bbox="630 706 1980 1079">For illustrative purposes only, suppose a patient is being fitted with the hearing aid system 602_B. During the fitting process, a piece of software may be executed on the mobile device 606 to interact with the patient wearing the hearing aid system 602_B. Such interaction includes sending audio test signals from the software executing on the mobile device 606 to the hearing aid system 602_B. Such software includes a user interface. The aural response from the patient is either automatically sent back to the testing software by the hearing aid system 602_B or is recorded manually into the testing software through the mobile device 606. Once sufficient aural responses are collected, in one embodiment, such aural responses are sent to a programming fitting server 616 through the cellular/mobile network 610 and the Internet 614; in another embodiment, such aural responses are formatted to form an audiogram before sending the information to the programming fitting server 616.</p> <p data-bbox="583 1096 848 1128">Eaton '458 at FIG. 6</p>

Claims	Prior Art
	 <p>FIG. 6</p> <p>A POSITA would have been motivated to combine the teachings of the DHA system with Anderson and Eaton to meet the claim limitations with a reasonable expectation of success for the reasons described above with respect to claim limitation [1b]. A POSITA would have understood Anderson '783's hearing test to teach receiving a user indication corresponding to the user's ability to hear the tones generated by the RPU because a POSITA would have understood such tests of hearing loss to rely on the user to indicate which audio tones are heard, for example, as disclosed by Eaton '458.</p> <p>To the extent that the DHA system, alone or in view of Anderson, does not disclose any elements of this claim limitation, it would have been obvious to a person having ordinary skill in the art to combine, with the teachings of DHA systems, those elements, as known to them or as disclosed by other prior art identified in the cover pleading.</p>
<p>4. A personal audio system as recited in claim 1, wherein the signal amplification is frequency-dependent amplification.</p>	<p>DHA systems, alone or in view of Anderson, expressly or inherently discloses “[a] personal audio system as recited in claim 1, wherein the signal amplification is frequency-dependent amplification.”</p> <p>For example, Anderson discloses that a program stored in the RPU “performs signal enhancement to compensate for a user’s hearing loss” using “appropriate gain vs. frequency parameters.” See, e.g., id. at</p>

Claims	Prior Art
	<p>27:13-17. A POSITA would have understood the disclosed “gain” to refer to a gain factor, indicating a level of signal amplification; and, accordingly, that the disclosed “gain vs. frequency parameters” for signal enhancement includes frequency-dependent gain factors used to perform frequency-dependent amplification. The following citations show examples of where Anderson ’783 discloses this limitation:</p> <p><i>See above at claim [1d].</i></p> <p>A POSITA would have been motivated to combine the teachings of the DHA system with Anderson to meet the claim limitations with a reasonable expectation of success for the reasons described above with respect to claim limitation [1b] and [1d].</p> <p>To the extent that the DHA system, alone or in view of Anderson, does not disclose any elements of this claim limitation, it would have been obvious to a person having ordinary skill in the art to combine, with the teachings of DHA systems, those elements, as known to them or as disclosed by other prior art identified in the cover pleading.</p>
<p>5[pre]. A personal audio system as recited in claim 1, wherein each of the earbuds includes a speaker, and</p>	<p>DHA systems, alone or in view of Anderson, expressly or inherently discloses “[a] personal audio system as recited in claim 1, wherein each of the earbuds includes a speaker.”</p> <p><i>See above at claim [1a].</i></p> <p>A POSITA would have been motivated to combine the teachings of the DHA system with Anderson to meet the claim limitations with a reasonable expectation of success for the reasons described above with respect to claim [1a].</p> <p>To the extent that the DHA system, alone or in view of Anderson, does not disclose any elements of this claim limitation, it would have been obvious to a person having ordinary skill in the art to combine, with the teachings of DHA systems, those elements, as known to them or as disclosed by other prior art identified in the cover pleading.</p>
<p>[5a] wherein the at least a portion of the customized audio output being generated by each of the ear buds is generated by the corresponding speaker.</p>	<p>DHA systems, alone or in view of Anderson, expressly or inherently discloses “wherein the at least a portion of the customized audio output being generated by each of the ear buds is generated by the corresponding speaker.”</p> <p><i>See above at claim [1g].</i></p> <p>A POSITA would have been motivated to combine the teachings of the DHA system with Anderson to meet the claim limitations with a reasonable expectation of success for the reasons described above with respect to claim [1a] and [1b].</p>

Claims	Prior Art
	<p>To the extent that the DHA system, alone or in view of Anderson, does not disclose any elements of this claim limitation, it would have been obvious to a person having ordinary skill in the art to combine, with the teachings of DHA systems, those elements, as known to them or as disclosed by other prior art identified in the cover pleading.</p>
<p>6. A personal audio system as recited in claim 1, wherein at least one of the earbuds in the pair of earbuds comprises an array of microphones configured for directional audio pickup.</p>	<p>DHA systems, alone or in view of Anderson, expressly or inherently discloses “[a] personal audio system as recited in claim 1, wherein at least one of the earbuds in the pair of earbuds comprises an array of microphones configured for directional audio pickup.”</p> <p>Phonak – Perseo Brochure at Page 23</p>  <p>Perseo BTEs – reliable beauty</p> <p>The Perseo BTE family includes a full choice of sizes, colors and power levels. All feature an elegant, award-winning design with a refined, ergonomic shape. They have a patented housing construction with bonded plastic and silicon to effectively keep out debris and moisture. Environmental protection extends further with patented, disposable microphone covers to ensure consistent, reliable performance, in all conditions. With dual-microphone systems this is particularly important since electronic matching cannot successfully compensate for performance reduction as a result of debris and moisture. A comprehensive range of colors for the housing also caters for specific tastes and trends.</p> <p>Original size: New small Perseo BTE 111/211</p>

Claims	Prior Art
	<p data-bbox="583 175 1627 212">Valente - Effect of type of noise and loudspeaker array on the performance at 399</p> <p data-bbox="638 228 1022 493"> <i>Recently, Phonak introduced its digital hearing aid in which one feature is an adaptive directional microphone. This adaptive microphone reportedly ensures identical sensitivity and phase characteristics of the dual microphones and a mechanism to <i>automatically</i> change the polar pattern from omnidirectional to any directional pattern while compensating for the changes in low-frequency response associated with each pattern. This hearing aid can be programmed</i> </p> <p data-bbox="583 516 1514 553">Phonak – The Hearing Review Article re Perseo (Sep. 5, 2003) at Page 2</p> <p data-bbox="638 570 1224 862"> <i>Perseo is a 20-channel digital hearing instrument. The default base-program is optimized for speech understanding in quiet environments. A second base-program with different frequency response settings and optional adaptive directional microphones provides additional support in noisy situations. The acoustic environment is continuously monitored and the appropriate base-program is automatically selected. This dual-base automatic function allows for independent fine tuning in one listening situation without compromising performance in the other. The automatic performance is designed to be easily personalized for the individual needs and preferences of the user.</i> </p> <p data-bbox="583 885 1176 922">Phonak – Annual Report (2003–2004) at p. 12</p> <p data-bbox="638 938 1085 1328"> <i>Within the FM business segment, the Group continued to leverage existing key technology into new applications. With SmartLink, Phonak has achieved another industry first by allowing hearing instrument users to access other communications systems. SmartLink combines a remote control, a transmitter, a directional microphone system and a Bluetooth link to mobile phone all in one single device. Phonak also presented EduLink, an innovative FM device to assist the learning development of children with Auditory Processing Disorder (APD). EduLink was perceived very positively within the medical community and recent classroom tests show favourable results.</i> </p>

Claims

Prior Art

Phonak – Collections 2002 Product List at Page 7



Adaptive digital AudioZoom (dAZ)



Fine-scale Noise Canceler (FNC)

Active and adaptive separation of speech from noise
The Speech-to-Noise Ratio (SNR) can only be optimised if the noise is spatially separated from the speech signal with directional microphone technology. As a pioneer in multi-microphone technologies, Phonak introduced AudioZoom, a feature which significantly improves listening and understanding in noisy situations, enhancing user satisfaction. Now, Phonak goes a step further in the refinement of multi-microphone technology with Adaptive digital AudioZoom: it seeks, finds, and suppresses noise automatically.

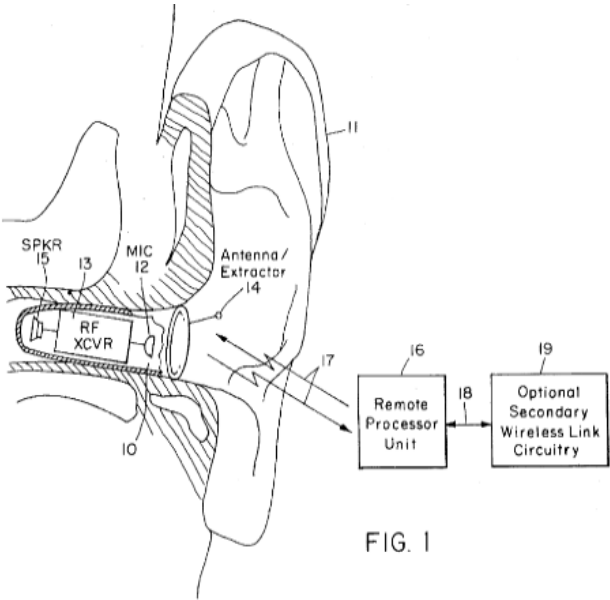
Fully automatic, fully adaptive
With its unique spatial processing system, Phonak's Adaptive digital AudioZoom can detect noise not only from directly behind, but anywhere in the near plane. Whether the hearing-impaired person is moving or the location of the noise source changes, the directional characteristic is constantly adapted to the listening situation. In this manner, the Adaptive digital AudioZoom automatically establishes the best possible SNR.

Twenty-band noise control
Claro constantly monitors the Speech-to-Noise Ratio (SNR) in all 20 bands. If a poor SNR is measured in certain bands, the Fine-scale Noise Canceler will reduce gain only in those bands. The result of this precision noise control is a selective and intelligent improvement of the SNR with only a marginal influence on the speech signal.

Precise and selective
All frequency bands are separately analysed and controlled. Gain is reduced selectively only in those bands where noise is present. This high-resolution analysis and gain control makes it possible to suppress noise with extremely high accuracy. Neighboring frequency bands without noise remain unaffected.

Intelligent
The magnitude of gain reduction is not constant. In frequency bands with higher SNRs and in bands which are important for speech intelligibility, gain is reduced only slightly. This intelligent control algorithm reduces the masking and distractive effects of background noise. Important speech information is preserved.

Alternatively, Anderson discloses this feature.

Claims	Prior Art
	<p data-bbox="583 175 894 207">Anderson '783 at FIG. 1</p>  <p data-bbox="995 748 1062 776">FIG. 1</p> <p data-bbox="583 846 945 878">Anderson '783 at 3:52–4:14</p> <p data-bbox="627 883 1986 1463">One preferred embodiment of the invention herein described is shown in block diagram form in FIG. 1. In this figure, the earpiece 10 shown worn in the ear 11 uses a standard completely-in-the-canal (CIC) housing, well known in the art, although many other earpiece housing types known in the art (e.g., behind-the-ear, or BTE) can be used. The earpiece shown in FIG. 1 comprises a microphone 12, RF transceiver 13 with an antenna that doubles as an earpiece extractor 14, and speaker 15. Although the dual-use antenna/extractor arrangement shown in FIG. 1 is often desirable, a separate antenna or antennae can be used in conjunction with many different types of earpiece extractors as appropriate for a particular application. Antennae and extractors can be disguised as ear hair or jewelry (e.g., earrings), or conductive filaments can be permanently implanted in the ear cartilage to act as antennae using a process similar to ear piercing. The earpiece 10 communicates with a remote processor unit (RPU) 16 via a primary two-way RF link 17, although many other wireless link media (e.g., ultrasonic or infrared) may be used instead. The RPU is typically worn under clothing (e.g., carried in a pocket or purse), but may also be worn in plain sight (e.g., on a belt) if desired. The RPU 16 may be connected (via wired or wireless means 18) to optional secondary wireless link circuitry 19 that allows wireless communication between the RPU and other sources of information (e.g., the general subscriber telephone network) via a secondary wireless link. Note that the optional secondary wireless link circuitry 19 may or may not be contained within the</p>

Claims	Prior Art
	<p>RPU case 70 of FIG. 7.</p> <p>Anderson '783 at 5:3–21</p> <p>During normal operation, speech signals from a nearby talker (and other signals in the ambient audio environment) are picked up by the earpiece wireless microphone 20 and transmitted to the RPU FM receiver 24. The RPU FM receiver 24 output level may be adjusted using the RPU FM receiver 24 volume control. The resulting electrical waveform representing signals from the nearby talker travels through the DPDT switch 25, which is set to the lower position as shown in FIG. 2, to the signal enhancer 26. The signal enhancer 26 may be, for example, a voice changer device that varies the pitch of a received speech signal according to the setting of pushbutton controls located on the RPU signal enhancer 26. The speech signal's pitch can then be raised or lowered as desired to help compensate for a user's hearing loss relative to a particular talker's voice characteristics. The modified speech signal travels from the RPU signal enhancer 26 to the RPU FM transmitter 27, and finally to the earpiece headset 21 receiver where the signal is converted to acoustic waves heard by the user.</p> <p>Anderson '783 at 6:8–25</p> <p>The user's voice (and other signals in the ambient audio environment) is picked up by the earpiece wireless microphone 20, transmitted to the RPU 23 FM receiver 24, sent through the DPDT switch 25 to the microphone input of the handset 28 (the handset's microphone transducer may be removed) and transmitted on the secondary wireless link 29 to, e.g., the general subscriber telephone network. Signals from the secondary wireless link 29 are received by the handset 28, sent from the speaker output of the handset 28 (the handset's speaker may be removed) via the DPDT switch 25 to the signal enhancer 26, FM transmitter 27 and finally to the earpiece headset 21 receiver where the signals are converted to acoustic waves heard by the user. Note that the RPU 23 may also include straight-forward connections to a "personal digital assistant" computer, or voice-operated devices connected to the FM receiver 24 output and not shown in FIG. 2.</p> <p>Anderson '783 at 22:11–61</p> <p>Stand-alone earpiece operation: If an earpiece cannot communicate with an RPU for any reason, stand-alone earpiece operation as a CIC hearing aid is generally desirable. Examples include occasions where the RPU must be turned off by pushing the RPU power-off key 78 (e.g., situations where operation of wireless transmitters is prohibited such as aircraft takeoffs and classified military briefings), RPU battery or other RPU electronics failure, misplacement of an RPU, and primary wireless link 17 interference. Note that a microphone located in the earpiece 10 is preferred to a body-worn microphone, since only an</p>

Claims	Prior Art
	<p>earpiece microphone 12 can be used to provide the desired default CIC hearing aid function (earpiece microphones also provide superior audio localization capabilities as the microphone moves with the user's head, and provide reduced wind and clothing noise compared with body-worn microphones). Parameters that control the stand-alone earpiece operation are stored in the earpiece CIC default parameter storage 870, which is accessed by the ASIC 845 to implement the desired stand-alone functions as described earlier. CIC default parameters include, but are not limited to, settings for the nominal and maximum allowable earpiece volume. Other information contained in the default parameter storage 870 includes, but is not limited to, the communication address and L/R (left/right) selection bit. To ensure that the CIC default parameter storage information cannot be changed inadvertently, changes are only allowed when the earpiece is not in use as, for example, when the earpiece battery 860 is being charged. A battery charging detector 840, that is sensitive to the high-level magnetic induction fields commonly used for earpiece battery charging, indicates to the ASIC 845 that the default parameter storage 870 values may be changed upon receipt of an appropriate command from the RPU. Specifically, to change information in the earpiece CIC default parameter storage 870, a user first activates the earpiece battery charging system, e.g., by placing the earpiece in a recharger/storage tray well known in the art. The RPU keyboard 946 is then used to enter a command to the DSP 948 that indicates the user's desire to change information in the default parameter storage 870, followed by the new information to be stored. The DSP 948 causes the RPU to transmit an interrogation with a reserved address code to which all earpieces being charged within the communication range of the RPU are receptive. The information that is to ultimately be loaded into the earpiece default parameter storage 870 is then transmitted by the RPU using the auxiliary bit field of many successive interrogations, said interrogations are received by the earpiece, and the desired information is loaded into the CIC default parameter storage 870 under control of the ASIC 845.</p> <p>Anderson '783 at 22:61–23:35</p> <p>Details of normal system operation: During normal operation of the hearing aid with wireless remote processor, sounds from the ambient environment are picked up by a microphone 12 in the earpiece 10 and transmitted with other information (e.g., data representing the earpiece battery 860 voltage as obtained by an A/D 865 converter and transmitted under control of the earpiece ASIC 845 in the auxiliary bit of many successive replies via the reply data register 825) over a primary two-way wireless link 17 to the RPU 16, where the audio signals are enhanced according to the user's needs by the RPU DSP 948. Signal enhancement is achieved via any of a number of techniques well known in the art that use the capabilities of the general-purpose RPU DSP 948. The enhanced audio signals may be combined with other information (e.g., a synthesized voice, generated within the DSP 948, warning that the aforementioned earpiece battery voltage is low) and transmitted from the RPU 16 over the primary wireless link 17 to the earpiece 10, where they are converted by a speaker 15 to sounds that can only be heard by the user 11. Art</p>

Claims	Prior Art
	<p>auxiliary speaker can be supplied in the RPU 16 as a peripheral 950 device that allows people other than the user to hear the same sounds as those produced by the earpiece speaker 15, if desired (as, e.g., when a person other than the user wishes to use the system's optional cellular telephone capabilities). A different, optional secondary two-way wireless link (shown as part of the RPU in FIG. 9 as secondary wireless link circuitry 944, or connected to the RPU by a wired or other, e.g. infrared, wireless link 18 and shown as optional secondary wireless link circuitry 19 in FIG. 1) can be used for communication between the RPU and a cellular telephone system or other sources of information. An RPU keyboard suitable for covert data entry 72, or voice recognition capabilities in the RPU (implemented, e.g., using any of a number of techniques well known in the art that use the general-purpose RPU DSP 948), can be used to control hearing aid parameters (such as amplification level set by the RPU DSP 948) and telephone dialing functions (e.g., automatic dialing of a predefined telephone number using information stored in the RPU DSP 948) in a manner that can be made imperceptible to a casual observer if desired.</p> <p>Anderson '783 at 23:36–42 Binaural, directional noise-canceling hearing aid with hearing protection: Two earpieces (one in each ear) and an RPU can be used in a binaural, directional noise-canceling hearing aid system with hearing protection features. Such a system provides some measure of protection for the residual hearing capability of a hearing-impaired user simultaneous with the hearing aid function.</p> <p>Anderson '783 at 25:29–30 The verbal commands are received as audio signals by the earpiece 10 and transmitted to the RPU 16.</p> <p>Anderson '783 at 26:54–65 Secondary link audio considerations: In some applications, a user may not wish to transmit signals from the ambient audio environment over the secondary wireless link. In such cases, an auxiliary earpiece microphone which favors the "bone conduction" path over the normal "air path" between the user's mouth and ear can be used. To save space, the auxiliary microphone may be housed in the same miniature case as the primary earpiece microphone. Alternatively, when two earpieces are used, unwanted ambient sources can be "nulled out" in the RPU DSP 948 by digital signal processing techniques well known in the art (e.g., the equivalent of speakerphone circuitry).</p> <p>Anderson '783 at 26:65–67 Note that such binaural signal processing techniques can also provide directivity and background noise</p>

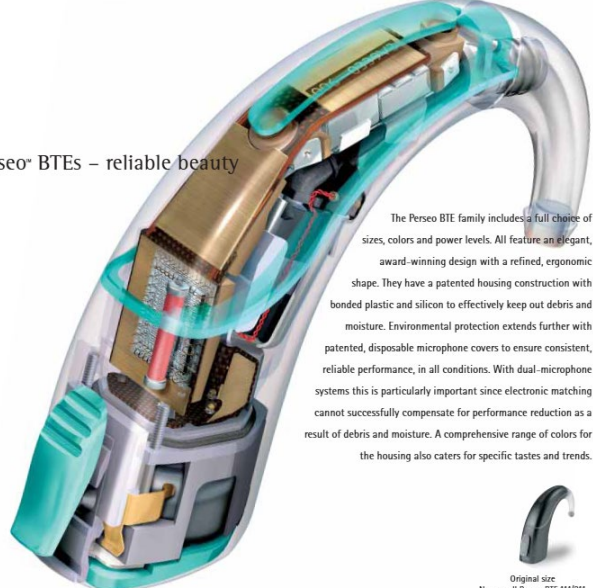
Claims	Prior Art
	<p>cancellation for users</p> <p>Anderson '783 at 27:1–3 having total hearing loss in one ear and partial loss in the other ear, as discussed earlier in connection with binaural hearing aid system operation.</p> <p>Alternatively, Jannard discloses this feature, for example, as explained in more detail in the corresponding portions of Exhibits B-1 and B-2. It would have been obvious to combine the DHA Systems with Anderson and Jannard for the same reasons set forth in Exhibits B-1 and B-2 with respect to combining Anderson and Jannard.</p> <p>A POSITA would have been motivated to combine the teachings of the DHA system with Anderson to meet the claim limitations with a reasonable expectation of success for the reasons described above with respect to claim [1a].</p> <p>To the extent that the DHA system, alone or in view of Anderson, does not disclose any elements of this claim limitation, it would have been obvious to a person having ordinary skill in the art to combine, with the teachings of DHA systems, those elements, as known to them or as disclosed by other prior art identified in the cover pleading.</p>
<p>7. A personal audio system as recited in claim 6, wherein direction of the directional audio pickup is in front of the user.</p>	<p>DHA systems, alone or in view of Anderson, expressly or inherently discloses “[a] personal audio system as recited in claim 6, wherein direction of the directional audio pickup is in front of the user.”</p> <p><i>See above at claim 6.</i></p> <p>A POSITA would have been motivated to combine the teachings of the DHA system with Anderson to meet the claim limitations with a reasonable expectation of success for the reasons described above with respect to claim [1a].</p> <p>To the extent that the DHA system, alone or in view of Anderson, does not disclose any elements of this claim limitation, it would have been obvious to a person having ordinary skill in the art to combine, with the teachings of DHA systems, those elements, as known to them or as disclosed by other prior art identified in the cover pleading.</p>
<p>8. A personal audio system as recited in claim 6, wherein direction of the directional audio pickup is outward from the user.</p>	<p>DHA systems, alone or in view of Anderson, expressly or inherently discloses “[a] personal audio system as recited in claim 6, wherein direction of the directional audio pickup is outward from the user.”</p> <p><i>See above at claim 6.</i></p>

Claims	Prior Art
	<p>A POSITA would have been motivated to combine the teachings of the DHA system with Anderson to meet the claim limitations with a reasonable expectation of success for the reasons described above with respect to claim [1a].</p> <p>To the extent that the DHA system, alone or in view of Anderson, does not disclose any elements of this claim limitation, it would have been obvious to a person having ordinary skill in the art to combine, with the teachings of DHA systems, those elements, as known to them or as disclosed by other prior art identified in the cover pleading.</p>
<p>9. A personal audio system as recited in claim 1, wherein at least one of the earbuds in the pair of earbuds comprises a plurality of microphones configured for audio pickup.</p>	<p>DHA systems, alone or in view of Anderson, expressly or inherently discloses “[a] personal audio system as recited in claim 1, wherein at least one of the earbuds in the pair of earbuds comprises a plurality of microphones configured for audio pickup.”</p> <p><i>See above at claim 6.</i></p> <p>A POSITA would have been motivated to combine the teachings of the DHA system with Anderson to meet the claim limitations with a reasonable expectation of success for the reasons described above with respect to claim [1a].</p> <p>To the extent that the DHA system, alone or in view of Anderson, does not disclose any elements of this claim limitation, it would have been obvious to a person having ordinary skill in the art to combine, with the teachings of DHA systems, those elements, as known to them or as disclosed by other prior art identified in the cover pleading.</p>
<p>10[pre]. A personal audio system as recited in claim 9, wherein the plurality of microphones includes at least:</p>	<p>DHA systems, alone or in view of Anderson, expressly or inherently discloses “[a] personal audio system as recited in claim 9, wherein the plurality of microphones includes at least.”</p> <p><i>See above at claim 6.</i></p> <p>A POSITA would have been motivated to combine the teachings of the DHA system with Anderson to meet the claim limitations with a reasonable expectation of success for the reasons described above with respect to claim [1a].</p> <p>To the extent that the DHA system, alone or in view of Anderson, does not disclose any elements of this claim limitation, it would have been obvious to a person having ordinary skill in the art to combine, with the teachings of DHA systems, those elements, as known to them or as disclosed by other prior art identified in the cover pleading.</p>

Claims	Prior Art
<p>[10a] a first microphone configured to receive audio pickup from the user; and</p>	<p>DHA systems, alone or in view of Anderson, expressly or inherently discloses “a first microphone configured to receive audio pickup from the user.”</p> <p><i>See above at claim 6.</i></p> <p>A POSITA would have been motivated to combine the teachings of the DHA system with Anderson to meet the claim limitations with a reasonable expectation of success for the reasons described above with respect to claim [1a].</p> <p>To the extent that the DHA system, alone or in view of Anderson, does not disclose any elements of this claim limitation, it would have been obvious to a person having ordinary skill in the art to combine, with the teachings of DHA systems, those elements, as known to them or as disclosed by other prior art identified in the cover pleading.</p>
<p>[10b] a second microphone configured to receive audio pickup outward away from the user.</p>	<p>DHA systems, alone or in view of Anderson, expressly or inherently discloses “a second microphone configured to receive audio pickup outward away from the user.”</p> <p><i>See above at claim 6.</i></p> <p>A POSITA would have been motivated to combine the teachings of the DHA system with Anderson to meet the claim limitations with a reasonable expectation of success for the reasons described above with respect to claim [1a].</p> <p>To the extent that the DHA system, alone or in view of Anderson, does not disclose any elements of this claim limitation, it would have been obvious to a person having ordinary skill in the art to combine, with the teachings of DHA systems, those elements, as known to them or as disclosed by other prior art identified in the cover pleading.</p>
<p>11. A personal audio system as recited in claim 10, wherein the audio pickup from at least the second microphone is used at least for noise cancelation.</p>	<p>DHA systems, alone or in view of Anderson and Nemirovski, expressly or inherently discloses “[a] personal audio system as recited in claim 10, wherein the audio pickup from at least the second microphone is used at least for noise cancelation.”</p>

Claims	Prior Art
	<p>Phonak – Perseo Brochure at Page 7</p> <div data-bbox="636 232 737 332" style="display: inline-block; width: 40px; height: 60px; background-color: #FFD700; margin-right: 10px;"></div> <div data-bbox="758 232 1205 297" style="display: inline-block; vertical-align: top;"> <h3>Convenience and Comfort</h3> <p>Perseo is easier to wear and easier to fit</p> </div> <div data-bbox="1230 232 1331 332" style="display: inline-block; width: 40px; height: 60px; background-color: #00AEEF; margin-right: 10px;"></div> <div data-bbox="1352 232 1816 297" style="display: inline-block; vertical-align: top;"> <h3>Clarity</h3> <p>Enhanced performance in multiple listening situations</p> </div> <div data-bbox="732 354 1822 748" style="display: flex; flex-wrap: wrap;"> <div style="width: 50%;"> <ul style="list-style-type: none"> ■ For true "wear and hear" convenience and comfort in all situations, Perseo automatic operation is managed by Phonak-exclusive PersonalLogic. ■ This breakthrough intelligence for automatic hearing instruments responds to acoustics, the listening situation and personal listening preferences. ■ Perseo is easily programmed for individual needs using the Listening Situation Manager. </div> <div style="width: 50%;"> <ul style="list-style-type: none"> ■ The Clarity Component Enhancement system realized in Digital Perception Processing² (DPP²) increases the clarity of speech and environmental sounds. ■ Rapid-tracking Adaptive digital AudioZoom (dAZ) is automatically activated to improve speech understanding in noisy situations. ■ For extra comfort and ease of listening, the Fine-scale Noise Canceler selectively reduces background noise in 20 channels. ■ Optional MicroLink FM provides the ultimate Speech-to-Noise Ratio enhancement. </div> </div> <p>Phonak – Perseo Brochure at Page 11</p> <ul style="list-style-type: none"> ■ With dual-base automatic the fitter has direct access to the two separate base-programs. Performance can be matched to each user's needs in quiet situations without disrupting performance in noisy situations and vice versa. The result is less fine-tuning effort than with single-base instruments. With single-base automatic time is wasted "juggling" a single set of fitting parameters while attempting to balance performance in multiple listening situations. ■ To further simplify fitting, PersonalLogic is pre-set according to typical needs and is easily re-programmed to individual preferences using the Listening Situation Manager (LSM). Ensuring that an automatic hearing system responds correctly in noisy situations where there is no speech is an important consideration for safety and to eliminate feelings of isolation. The LSM provides this functionality with a single mouse click. Matching the system's sensitivity to noise with that of the client's, and controlling the auditory contrast between quiet and noisy settings is equally easy.

Claims	Prior Art
	<p data-bbox="583 175 1062 207">Phonak – Perseo Brochure at Page 13</p> <div data-bbox="627 264 1482 308">  <h2 data-bbox="636 264 1474 308">Enhanced sound quality and Speech-to-Noise Ratio</h2> </div> <p data-bbox="726 362 1062 383">Clarity Component Enhancement – DPP²</p> <p data-bbox="726 401 1719 615">Speech clarity depends on both its acoustic properties and the listener's auditory perception. Conventional spectral enhancement systems consider only the acoustics of sound, not how it will be perceived by the listener. Perseo goes well beyond this with DPP², the new generation of Digital Perception Processing. DPP² with Clarity Component Enhancement is individualized for each wearer. It analyzes the acoustic and the psychoacoustic characteristics of incoming sound in 20 critical bands to identify accurately the components that contribute to clarity. It then selectively amplifies them to optimize clarity for the wearer.</p> <p data-bbox="583 638 1024 670">Phonak – Perseo Brochure at p. 14</p> <p data-bbox="636 688 840 716">Rapid-tracking dAZ</p> <p data-bbox="636 742 1056 1122">Phonak pioneered the development of directional microphone systems. Not surprisingly, Perseo dAZ with new soft switching brings even better performance. Perseo dAZ is automatically activated in appropriate situations and immediately reduces noise from beside and behind the wearer. The wearer experiences an immediate reduction of background noise and enhanced speech clarity. Rapid tracking ensures that even fast moving noise sources are controlled.</p>

Claims	Prior Art
	<p data-bbox="583 175 1062 207">Phonak – Perseo Brochure at Page 23</p> <div data-bbox="625 224 1260 824">  <p data-bbox="636 383 905 407">Perseo® BTEs – reliable beauty</p> <p data-bbox="968 443 1255 695">The Perseo BTE family includes a full choice of sizes, colors and power levels. All feature an elegant, award-winning design with a refined, ergonomic shape. They have a patented housing construction with bonded plastic and silicon to effectively keep out debris and moisture. Environmental protection extends further with patented, disposable microphone covers to ensure consistent, reliable performance, in all conditions. With dual-microphone systems this is particularly important since electronic matching cannot successfully compensate for performance reduction as a result of debris and moisture. A comprehensive range of colors for the housing also caters for specific tastes and trends.</p> <p data-bbox="1140 800 1247 824">Original size New small Perseo BTE 111/211</p> </div> <p data-bbox="583 846 1514 878">Phonak – The Hearing Review Article re Perseo (Sep. 5, 2003) at Page 2</p> <p data-bbox="636 898 1224 1190">Perseo is a 20-channel digital hearing instrument. The default base-program is optimized for speech understanding in quiet environments. A second base-program with different frequency response settings and optional adaptive directional microphones provides additional support in noisy situations. The acoustic environment is continuously monitored and the appropriate base-program is automatically selected. This dual-base automatic function allows for independent fine tuning in one listening situation without compromising performance in the other. The automatic performance is designed to be easily personalized for the individual needs and preferences of the user.</p> <p data-bbox="583 1214 1119 1247">Phonak – Microlink User Guide at Page 2</p> <p data-bbox="636 1266 1161 1425">The desired signal – for example, a speaker’s voice or the sound from a TV – is picked up by the microphone which is placed near to the source of the sound. The transmitter sends an FM signal to the receiver, which immediately relays this to the hearing instrument for you to hear – crystal clear!</p>

Claims

Prior Art

Phonak – Collections 2002 Product List at Page 7



Adaptive digital AudioZoom (dAZ)



Fine-scale Noise Canceler (FNC)

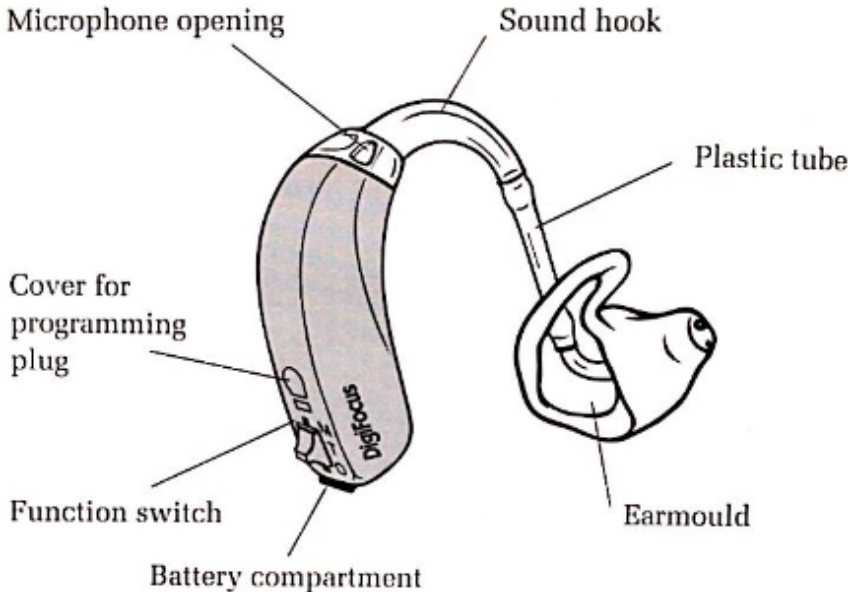
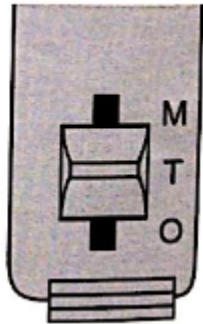
Active and adaptive separation of speech from noise
The Speech-to-Noise Ratio (SNR) can only be optimised if the noise is spatially separated from the speech signal with directional microphone technology. As a pioneer in multi-microphone technologies, Phonak introduced AudioZoom, a feature which significantly improves listening and understanding in noisy situations, enhancing user satisfaction. Now, Phonak goes a step further in the refinement of multi-microphone technology with Adaptive digital AudioZoom: it seeks, finds, and suppresses noise automatically.

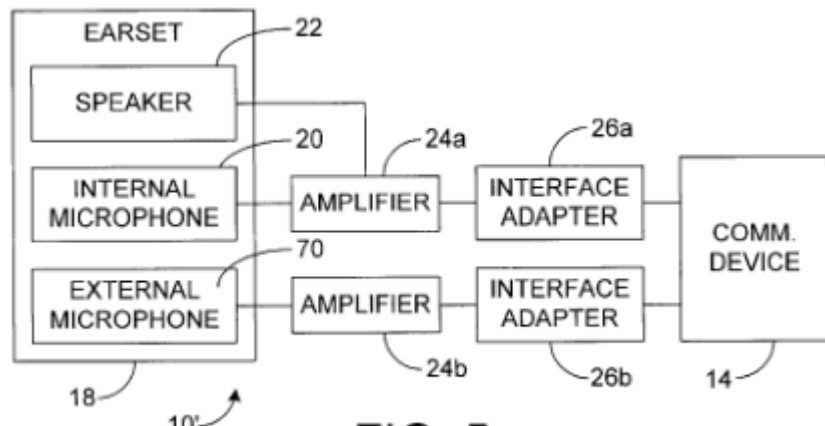
Fully automatic, fully adaptive
With its unique spatial processing system, Phonak's Adaptive digital AudioZoom can detect noise not only from directly behind, but anywhere in the near plane. Whether the hearing-impaired person is moving or the location of the noise source changes, the directional characteristic is constantly adapted to the listening situation. In this manner, the Adaptive digital AudioZoom automatically establishes the best possible SNR.

Twenty-band noise control
Claro constantly monitors the Speech-to-Noise Ratio (SNR) in all 20 bands. If a poor SNR is measured in certain bands, the Fine-scale Noise Canceler will reduce gain only in those bands. The result of this precision noise control is a selective and intelligent improvement of the SNR with only a marginal influence on the speech signal.

Precise and selective
All frequency bands are separately analysed and controlled. Gain is reduced selectively only in those bands where noise is present. This high-resolution analysis and gain control makes it possible to suppress noise with extremely high accuracy. Neighboring frequency bands without noise remain unaffected.

Intelligent
The magnitude of gain reduction is not constant. In frequency bands with higher SNRs and in bands which are important for speech intelligibility, gain is reduced only slightly. This intelligent control algorithm reduces the masking and distractive effects of background noise. Important speech information is preserved.

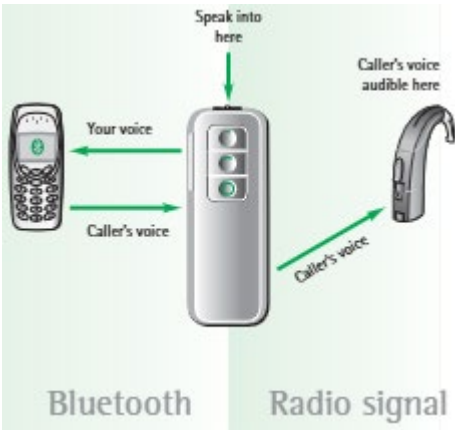
Claims	Prior Art
	<p data-bbox="583 175 1446 207">Oticon – DigiFocus Compact User Instructional Brochure at Page 4</p> <div data-bbox="632 228 1474 816">  </div> <p data-bbox="583 846 1446 878">Oticon – DigiFocus Compact User Instructional Brochure at Page 4</p> <p data-bbox="638 898 856 930">Function switch</p> <p data-bbox="638 941 1262 1008">M: Hearing instrument switched on and microphone connected. Normal use.</p> <p data-bbox="638 1049 1308 1263">T: Hearing instrument switched on and telecoil connected. Where a telemagnetic loop system has been installed or when you talk on the phone (in many countries), you may then be able to use your hearing instruments with minimum disturbance from background noise.</p> <p data-bbox="638 1300 1115 1333">O: Hearing instrument switched off.</p> <div data-bbox="1367 987 1568 1308">  </div> <p data-bbox="583 1365 974 1398">Anderson '783 at 23:65–24:14</p> <p data-bbox="632 1403 1940 1469">Background noise and competing talker cancellation are achieved using any of a number of techniques and circuits well known in the art, e.g., variations of the Least Mean Squares (LMS) technique as</p>

Claims	Prior Art
	<p>implemented in the Motorola DSP56200 LMS filter integrated circuit, or the equivalent function implemented in the general-purpose RPU DSP 948. Note that for all binaural signal processing functions, data contained in replies from the left earpiece (i.e., audio data from the left earpiece microphone) is processed in the RPU DSP 948 by combining said data with other data contained in replies from the right earpiece (i.e., audio data from the right earpiece microphone). The processed audio data transmitted from the RPU to the left and right earpiece speakers via RPU interrogations is therefore a function of data received from both the left and right earpiece microphones via earpiece replies.</p> <p>Nemirovski '379 at FIG. 5</p>  <p style="text-align: center;">FIG. 5</p> <p>Nemirovski '379 at 2:59–60</p> <p>FIG. 5 is a block diagram of an ear microphone assembly according to another aspect of the invention; and</p> <p>Nemirovski '379 at 3:6–23</p> <p>The present invention is directed to a communication system, including a communications device and an ear microphone, or earset, assembly. The communications system allows a user to speak with a remotely located person. The ear microphone assembly, also referred to herein as an earset, includes a microphone and a speaker supported by a housing. The housing is retained by the ear of the user and allows for hands free use of the communication system while carrying on a conversation with the remotely located person. The microphone is directed towards the ear canal of the user and detects sounds emanating or coming out of the ear (or air pressure changes occurring within the ear) to accurately and reliably detect the speech of</p>

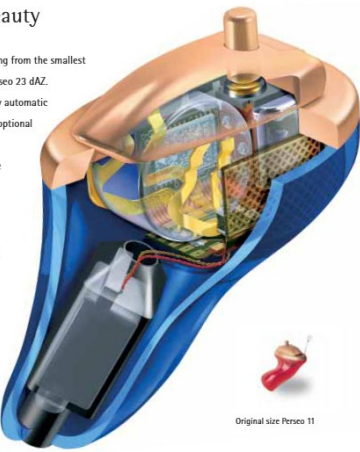
Claims	Prior Art
	<p>the user. The position of the microphone may also help to reduce detection of externally generated sounds. The ear microphone assembly allows for separation of a speech input device for a cellular, or wireless, telephone from an RF transceiver of the telephone.</p> <p>Nemirovski '379 at 4:27–41</p> <p>The communications system 10 includes an ear microphone assembly, or earset 18, having a microphone 20 and a speaker 22 supported by a housing 30 (FIG. 3). The physical arrangement of the earset 18 will be described in more detail below. The microphone 20 is used to detect sounds emanating or coming from the ear (sometimes referred to as air pressure changes within the ear) of the user that result from speech of the user. The microphone 20 converts those detections into an electrical signal. The electrical signal is optionally amplified by an amplifier 24. Next, the amplified electrical signal is received by the communication device 14 via an appropriate interface adapter 26. The interface adapter 26 can be, for example, a jack or other electrical connector received by a corresponding receptacle of the communication device 14.</p> <p>Nemirovski '379 at 9:33–39</p> <p>The external microphone 70 is disposed on a exterior portion 72 of the housing 30 and is pointed away from the wearer of the earset 18 assembly. The external microphone 70 can be secured to the housing 30 with adhesive and can be disposed in a recess defined by the exterior portion 72. Similar to the internal microphone 20, the external microphone 70 is small (e.g., less than two mm in size).</p> <p>Nemirovski '379 at 9:45–60</p> <p>FIG. 5 is a block diagram that illustrates a system 10' for detecting speech of the user. The system 10' includes the earset 18 having the internal microphone 20, the speaker 22 and amplifier 24a as described in greater detail above. In addition, the earset 18 of the system 10' has the external microphone 70. As one skilled in the art will appreciate, the use of the external microphone 70 is optional. The external microphone 70 is disposed on the housing 30 (FIG. 3) and generally points away from the person wearing the earset 18. As a result, the external microphone 70 detects air pressure changes located outside the ear. These air pressure changes are mainly the result of externally generated noise (other people's speech, environmental noise, etc.) and, in part, to the speech of the person wearing the earset 18, when the speech travels from the mouth of the user and through the air to the external microphone 70.</p>

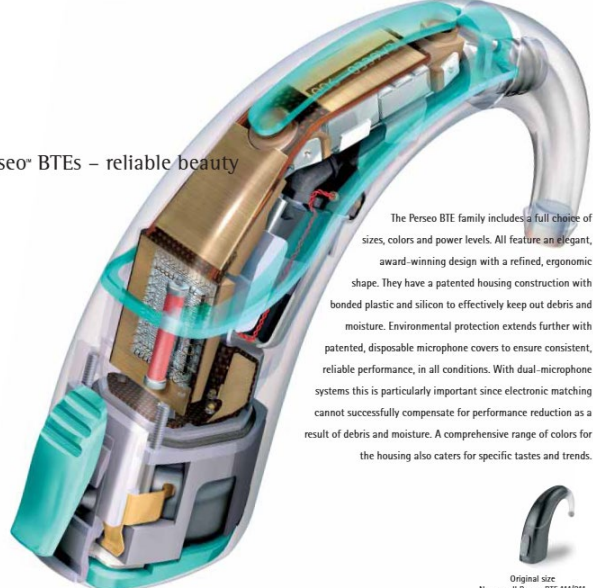
Claims	Prior Art
	<p>Nemirovski '379 at 10:8–13</p> <p>The unique signals generated by the internal and external microphone 20, 70 can be used to distinguish between user produced speech and externally produced sounds. In alternative configurations of the system 10', the system 10' can be adapted to cancel externally generated noise from the detections of the internal microphone 20.</p> <p>Nemirovski '379 at 10:14–26</p> <p>When the user speaks both the internal microphone 20 and the external microphone 70 detects the speech. Sounds generated by a source other than the user (e.g., another person speaking, a passing vehicle, machinery, the hum of an air conditioner or computer, etc.) will be referred to herein as noise. Both the external microphone 70 and the internal microphone 20 detect the noise. Accordingly, the communication device 14 (or other processing circuit) can be programmed to isolate the user's speech by identifying (i.e., separating) the speech component and the noise component. As a result, a relatively noise-free speech signal can be produced which can be transmitted to the communication network 12 by the communication device 14.</p> <p>Nemirovski '379 at 10:27–47</p> <p>In one embodiment, signal separation can be made by solving equations 1 and 2, identified below, for a sound component signal (“S”) and a noise component signal (“N”). An impulse response of the internal microphone (F_{MI}) and an impulse response of the external microphone (F_{ME}) reflect the sensitivity of the microphones, the frequency response of the microphones, resonating characteristics of the ear and housing, etc. Y_I represents the overall signal produced by the internal microphone and Y_E represents the overall signal produced by the external microphone. Assuming that the entire system is linear (the system features are additive and homogeneous), the system can be represented by equations 1 and 2 as follows (the asterisk, or “*”, represents a convolution operator):</p> <p>[Equation]</p> $F_{ME} * S + F_{ME} * N = Y_E \quad \text{Eq. 2}$ <p>By using the output signals Y_I and Y_E and the known impulse responses F_{MI} and F_{ME}, the two equations can be solved in real time for the two unknown signals S and N.</p> <p>As disclosed in the above–cited citations, Nemirovski '379 discloses a first microphone for picking up audio from a user (the internal microphone), a second microphone for picking up audio generated outward away from the user (the external microphone), and that the audio is used for noise cancelation. A POSITA would thus have understood that the audio picked up by Nemirovski '379's second (external) microphone is used for noise cancelation. A POSITA would have been motivated to combine the DHA system and Nemirovski '379</p>

Claims	Prior Art
	<p>to meet the claim limitations with a reasonable expectation of success for at least the reasons described above in reference to claim [1a]. The DHA system discloses that its hearing aid system can performs noise cancellation to enhance speech and Nemirovski '379 discloses one such example of a known technique for producing “a relatively noise-free speech signal” that includes “cancel[ing] externally generated noise from the detections of the internal microphone.” See id. at 10:11-14, 24-25. It therefore would have been obvious to a POSITA to employ Nemirovski '379’s noise cancellation techniques in the DHA system. The features of Nemirovski 379’s earset would have been readily combinable with the Phonak hearing aid by a POSITA with a reasonable expectation of success. Both systems disclose ear-worn devices for use with a communication system and can also include noise canceling functionality. Both devices also share other analogous and familiar computer components, such as microphones, speakers, and amplifiers. Accordingly, combining their respective teachings would involve combining known components to yield predictable results.</p> <p>To the extent that the DHA system, alone or in view of Anderson and Nemirovski, does not disclose any elements of this claim limitation, it would have been obvious to a person having ordinary skill in the art to combine, with the teachings of DHA systems, those elements, as known to them or as disclosed by other prior art identified in the cover pleading.</p>
<p>12. A personal audio system as recited in claim 11, wherein at least a portion of the audio pickup from the second microphone is subtracted from the audio pickup from the first microphone.</p>	<p>DHA systems, alone or in view of Anderson and Nemirovski, expressly or inherently discloses “[a] personal audio system as recited in claim 11, wherein at least a portion of the audio pickup from the second microphone is subtracted from the audio pickup from the first microphone.”</p> <p><i>See above at claim 11.</i></p> <p>A POSITA would have been motivated to combine the teachings of the DHA system with Anderson and Nemirovski to meet the claim limitations with a reasonable expectation of success for the reasons described above with respect to claims [1a] and 11.</p> <p>To the extent that the DHA system, alone or in view of Anderson and Nemirovski, does not disclose any elements of this claim limitation, it would have been obvious to a person having ordinary skill in the art to combine, with the teachings of DHA systems, those elements, as known to them or as disclosed by other prior art identified in the cover pleading.</p>
<p>13. A personal audio system as recited in claim 12, wherein each of the ear buds includes at least a speaker and is configured to generate at least a portion of the</p>	<p>DHA systems, alone or in view of Anderson, expressly or inherently discloses “[a] personal audio system as recited in claim 12, wherein each of the ear buds includes at least a speaker and is configured to generate at least a portion of the customized audio output for the user via the speaker.”</p>

Claims	Prior Art
<p>customized audio output for the user via the speaker.</p>	<p>Phonak - SmartLink SX User Guide at p. 32</p>  <p>Phonak - SmartLink SX User Guide at p. 39</p> <p><u>Answering a call</u></p> <ul style="list-style-type: none"> • Speak directly into the SmartLink. You hear the caller's voice in both hearing instruments. <p>Phonak - SmartLink SX User Guide at p. 58</p> <p>The user benefit of FM is an improved signal-to-noise ratio, which leads to an optimized speech understanding even in very difficult listening situations.</p>

Claims	Prior Art
	<p data-bbox="583 175 1010 207">Phonak – Perseo Brochure at p. 8</p>  <p data-bbox="583 846 1062 878">Phonak – Perseo Brochure at Page 13</p> <h3 data-bbox="636 935 1478 976">Enhanced sound quality and Speech-to-Noise Ratio</h3> <p data-bbox="726 1029 1062 1053">Clarity Component Enhancement – DPP²</p> <p data-bbox="726 1068 1719 1284">Speech clarity depends on both its acoustic properties and the listener's auditory perception. Conventional spectral enhancement systems consider only the acoustics of sound, not how it will be perceived by the listener. Perseo goes well beyond this with DPP², the new generation of Digital Perception Processing. DPP² with Clarity Component Enhancement is individualized for each wearer. It analyzes the acoustic and the psychoacoustic characteristics of incoming sound in 20 critical bands to identify accurately the components that contribute to clarity. It then selectively amplifies them to optimize clarity for the wearer.</p>

Claims	Prior Art
	<p>Phonak – Perseo Brochure at Page 22</p> <h1 data-bbox="655 224 1255 376">The beauty of high technology</h1> <p data-bbox="703 376 961 402">Perseo[®] ITEs – hidden beauty</p> <div data-bbox="630 487 655 535">2</div> <p data-bbox="703 425 1003 587">The Perseo family includes a full line of ITEs ranging from the smallest Perseo 11 CIC to the elegant dual-microphone Perseo 23 dAZ. All CIC and ITE instruments, like the BTEs, are fully automatic and have PersonalLogic. A TacTronic switch is optional on all ITE models and easily added after fitting by exchanging battery covers. This is a very simple and practical means of providing extra control and additional programs.</p> <p data-bbox="703 613 919 695">In many countries, ITEs are available with Phonak eShell. These electronically created shells are hypoallergenic, exceptionally durable and more easily retained in the wearer's ear.</p>  <p data-bbox="1150 787 1234 799">Original size Perseo 11</p>

Claims	Prior Art
	<p data-bbox="583 175 1062 207">Phonak – Perseo Brochure at Page 23</p> <div data-bbox="625 224 1260 824">  <p data-bbox="636 381 905 406">Perseo® BTEs – reliable beauty</p> <p data-bbox="968 440 1255 691">The Perseo BTE family includes a full choice of sizes, colors and power levels. All feature an elegant, award-winning design with a refined, ergonomic shape. They have a patented housing construction with bonded plastic and silicon to effectively keep out debris and moisture. Environmental protection extends further with patented, disposable microphone covers to ensure consistent, reliable performance, in all conditions. With dual-microphone systems this is particularly important since electronic matching cannot successfully compensate for performance reduction as a result of debris and moisture. A comprehensive range of colors for the housing also caters for specific tastes and trends.</p> <p data-bbox="1140 800 1247 821">Original size New small Perseo BTE 111/211</p> </div>

Claims

Prior Art

Phonak – Perseo Brochure at Page 26

Perseo 11



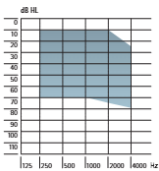
||| Digital microphone

OPTIONS

- ||| Tactronic switch
- ||| Control beeps

Ear Simulator data
45 dB max. gain
119 dB SPL max. output

2 cc data
35 dB max. gain
108 dB SPL max. output



Perseo 12



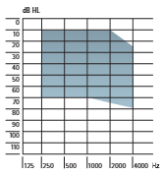
||| Digital microphone
||| Volume control

OPTIONS

- ||| Tactronic switch
- ||| Control beeps

Ear Simulator data
50 dB max. gain
119 dB SPL max. output

2 cc data
40 dB max. gain
108 dB SPL max. output



Perseo 22



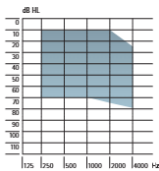
||| Telecoil

OPTIONS

- ||| Tactronic switch
- ||| Control beeps
- ||| Remote control WatchPlot2 and SoundPlot

Ear Simulator data
58 dB max. gain
118 dB SPL max. output

2 cc data
48 dB max. gain
108 dB SPL max. output



Claims

Prior Art

Phonak – Perseo Brochure at Page 27

Perseo 23 dAZ



■ Adaptive digital AudioZoom

■ Telecoil

■ Tactronic switch

■ Remote control WatchPlot2 and SoundPilot

■ Control beeps

■ Ear Simulator data

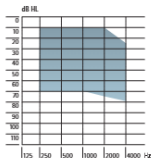
58 dB max. gain

118 dB SPL max. output

2 cc data

45 dB max. gain

108 dB SPL max. output



Perseo 111 dAZ
Perseo 111 dAZ FM



■ FM solutions

■ Adaptive digital AudioZoom

■ Telecoil

■ Tactronic switch

■ Remote control WatchPlot2 and SoundPilot

■ Control beeps

■ Ear Simulator data

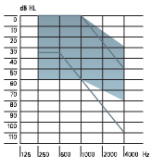
54 dB max. gain

127 dB SPL max. output

2 cc data

47 dB max. gain

121 dB SPL max. output



Perseo 211 dAZ
Perseo 211 dAZ FM



■ FM solutions

■ Adaptive digital AudioZoom

■ Telecoil

■ Tactronic switch

■ Remote control WatchPlot2 and SoundPilot

■ Control beeps

■ Ear Simulator data

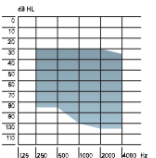
62 dB max. gain

130 dB SPL max. output

2 cc data

55 dB max. gain

125 dB SPL max. output



Technical Data

Claims

Prior Art

Phonak – Perseo Brochure at Page 28

Perseo 311 dAZ Forte
Perseo 311 dAZ Forte FM



Remote controls



FM Transmitters



- FM solutions
- Adaptive digital AudioZoom
- Telecoil
- Tactronic Switch

OPTIONS

- Control beeps
- Remote control WatchPilot2 and SoundPilot

Ear Simulator data

72 dB max. gain
135 dB SPL max. output

2 cc data*

85 dB max. gain
130 dB SPL max. output

WatchPilot™2



HandyMic™





SoundPilot™



TelCom™



Campus S™



Valente - Effect of type of noise and loudspeaker array on the performance at 398

Differences in performance between omnidirectional and directional microphones were evaluated between two loudspeaker conditions (single loudspeaker at 180°; diffuse using eight loudspeakers set 45° apart) and two types of noise (steady-state HINT noise; R-Space™ restaurant noise). Twenty-five participants were fit bilaterally with Phonak Perseo hearing aids using the manufacturer's recommended procedure. After wearing the hearing aids for one week, the

Phonak – Microlink User Guide at Page 2

The desired signal – for example, a speaker's voice or the sound from a TV – is picked up by the microphone which is placed near to the source of the sound. The transmitter sends an FM signal to the receiver, which immediately relays this to the hearing instrument for you to hear – crystal clear!

Claims	Prior Art																																																																																
	<p>Phonak – Claro 211 Product Information at Page 2</p> <p>Principle of FM operation</p> <p>The transmitter unit is worn by the person speaking, placed in front of the sound source or even connected directly to TV or HiFi. The transmitter unit collects the desired signal and transmits it via radio waves directly to the FM module integrated in the hearing instrument. The incoming FM signal received by the hearing instrument is processed either in accordance with the audiological setting of the instrument, providing consistent sound quality or according to an individual FM program.</p> <p>Phonak – Collections 2002 Product List at Page 5</p> <div><p>Elite</p><ul style="list-style-type: none">- Claro's digital sound processing with a human perspective (Digital Perception Processing)- AutoSelect: Excellent speech intelligibility in any listening situation- Adaptive digital AudioZoom- Fine-scale Noise Canceler: High-resolution 20-band noise control<table><tr><th>Digital</th><th></th><th></th><th></th><th></th><th></th><th></th><th></th><th></th><th></th></tr><tr><th>ITE</th><th></th><th></th><th></th><th>BTE</th><th></th><th></th><th></th><th>Remote controls</th><th></th></tr><tr><td>Claro 11</td><td></td><td></td><td></td><td>Claro 111 dAZ</td><td></td><td></td><td></td><td>SoundPilot FM</td><td></td></tr><tr><td>Claro 21 dAZ</td><td></td><td></td><td></td><td>Claro 211 dAZ</td><td></td><td></td><td></td><td>WatchPilot Elegance</td><td></td></tr><tr><td>Claro 22</td><td></td><td></td><td></td><td>Claro 311 dAZ</td><td></td><td></td><td></td><td>WatchPilot Classic</td><td></td></tr><tr><td></td><td></td><td></td><td></td><td>Claro 111 dAZ FM</td><td></td><td></td><td></td><td></td><td></td></tr><tr><td></td><td></td><td></td><td></td><td>Claro 211 dAZ FM</td><td></td><td></td><td></td><td></td><td></td></tr><tr><td></td><td></td><td></td><td></td><td>Claro 311 dAZ FM</td><td></td><td></td><td></td><td></td><td></td></tr></table><p>Elite Line features:</p><ul style="list-style-type: none">- Clear, natural sound with precise, fully automatic volume control- Significantly improved speech intelligibility even in difficult listening situations- Design innovations for all system components- High level of user convenience with automatic hearing program selection- Innovative high-tech solutions for Claro system control- State-of-the-art hearing computers for the optimised solution to a wide range of hearing problems- FM receiver fully integrated in the BTE case (Claro dAZ FM)</div>	Digital										ITE				BTE				Remote controls		Claro 11				Claro 111 dAZ				SoundPilot FM		Claro 21 dAZ				Claro 211 dAZ				WatchPilot Elegance		Claro 22				Claro 311 dAZ				WatchPilot Classic						Claro 111 dAZ FM										Claro 211 dAZ FM										Claro 311 dAZ FM					
Digital																																																																																	
ITE				BTE				Remote controls																																																																									
Claro 11				Claro 111 dAZ				SoundPilot FM																																																																									
Claro 21 dAZ				Claro 211 dAZ				WatchPilot Elegance																																																																									
Claro 22				Claro 311 dAZ				WatchPilot Classic																																																																									
				Claro 111 dAZ FM																																																																													
				Claro 211 dAZ FM																																																																													
				Claro 311 dAZ FM																																																																													

Claims

Prior Art

Phonak – Collections 2002 Product List at Page 6

Elite


Digital Perception Processing (DPP)


AutoSelect (AS)

New Phonak sound:
unmatched clarity
Digital Perception Processing (DPP) is based on a cochlear sound analysis across 20 critical bands – exactly like the analysis performed by the healthy human ear. But Clari gets a step further than previous overlapping multi-band filter systems. DPP models the active, load-ress-dependent cochlear signal processing in all 20 critical bands.

20 critical bands:
overlapping and dynamic
The cochlear filters are overlapping, but the degree of overlap is dependent on the frequency, level, and bandwidth of the sound signal. This is precisely how Clari's filters work. Digital Perception Processing computes the cochlear excitation patterns and thus determines sound perception. DPP is the first to mimic the healthy human ear in this way.

Automatic, intelligent program selection
Clari listens and automatically activates the optimal DPP algorithm for the acoustic environment: Program 1 with QuietAdapt for quiet situations or Program 2 with NoiseAdapt and noise control for noisy environments.

In quiet situations
QuietAdapt DPP is designed for clear hearing in quiet surroundings. QuietAdapt controls gain in 20 critical bands to correct for loudness and provide a natural sound quality. By preserving the perception patterns of complex sounds, QuietAdapt assures that low-level signal components are not over-

amplified to the point where they would adversely affect speech intelligibility. This allows Clari to deliver only those sound components that would be naturally audible.

In noise
NoiseAdapt DPP is the ideal program for maximizing audibility in noisy environments. The special NoiseAdapt algorithm works with the Fine-scale Noise Canceller – and if available with the Adaptive digital AudioZoom feature – to optimize speech understanding in noise.



Claims

Prior Art

Phonak – Collections 2002 Product List at Page 7



Adaptive digital AudioZoom (dAZ)



Fine-scale Noise Canceler (FNC)

Active and adaptive separation of speech from noise
The Speech-to-Noise Ratio (SNR) can only be optimised if the noise is spatially separated from the speech signal with directional microphone technology. As a pioneer in multi-microphone technologies, Phonak introduced AudioZoom, a feature which significantly improves listening and understanding in noisy situations, enhancing user satisfaction. Now, Phonak goes a step further in the refinement of multi-microphone technology with Adaptive digital AudioZoom: it seeks, finds, and suppresses noise automatically.

Fully automatic, fully adaptive
With its unique spatial processing system, Phonak's Adaptive digital AudioZoom can detect noise not only from directly behind, but anywhere in the near plane. Whether the hearing-impaired person is moving or the location of the noise source changes, the directional characteristic is constantly adapted to the listening situation. In this manner, the Adaptive digital AudioZoom automatically establishes the best possible SNR.

Twenty-band noise control
Claro constantly monitors the Speech-to-Noise Ratio (SNR) in all 20 bands. If a poor SNR is measured in certain bands, the Fine-scale Noise Canceler will reduce gain only in those bands. The result of this precision noise control is a selective and intelligent improvement of the SNR with only a marginal influence on the speech signal.

Precise and selective
All frequency bands are separately analysed and controlled. Gain is reduced selectively only in those bands where noise is present. This high-resolution analysis and gain control makes it possible to suppress noise with extremely high accuracy. Neighboring frequency bands without noise remain unaffected.

Intelligent
The magnitude of gain reduction is not constant. In frequency bands with higher SNRs and in bands which are important for speech intelligibility, gain is reduced only slightly. This intelligent control algorithm reduces the masking and distractive effects of background noise. Important speech information is preserved.

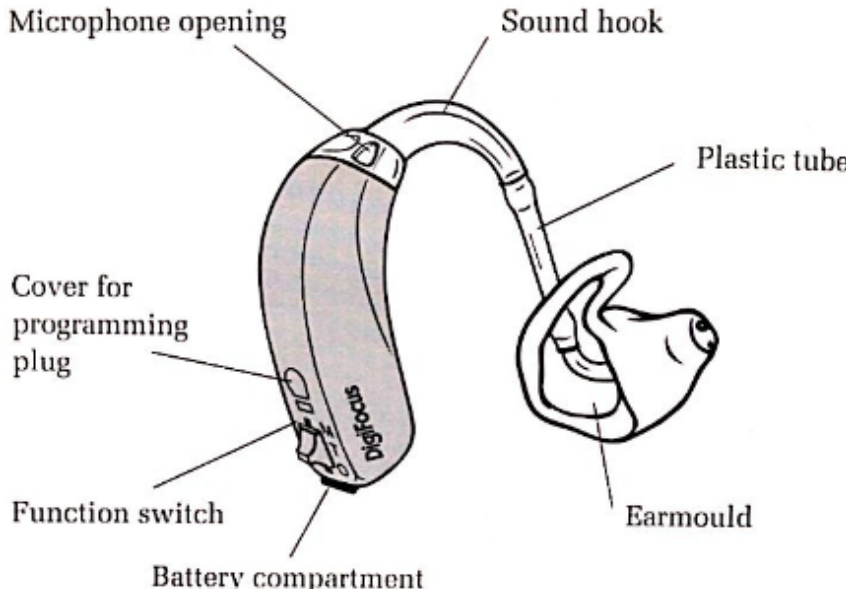
Claims

Prior Art

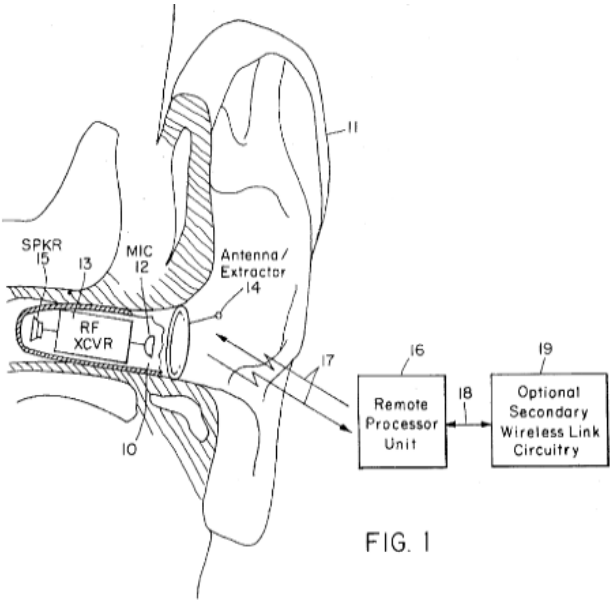
Phonak – Collections 2002 Product List at Page 8

	Claro 11 (CIC)	Claro 21 dAZ	Claro 22
ITE instruments			
			
Peak values G max./OSPL 90	45 dB/117 dB SPL	58 dB/117 dB SPL	58 dB/117 dB SPL
Features	World's first digital microphone TacTronic switch standard	In-the-Ear hearing computer with Adaptive digital Audio/oom (dAZ)	Semi-modular Claro In-the-Canal instrument
Description	Digital Perception Processing (DPP) AutoSelect (AS) Fine-scale Noise Canceler (FNC) 3 hearing programs + "mute" Battery size 10 TacTronic program selection switch	Digital Perception Processing (DPP) AutoSelect (AS) Fine-scale Noise Canceler (FNC) Adaptive digital Audio/oom (dAZ) 3 hearing programs + "mute" Battery size 312 Telecoil	Digital Perception Processing (DPP) AutoSelect (AS) Fine-scale Noise Canceler (FNC) 3 hearing programs + "mute" Battery size 312 Telecoil
Options	Rotary volume control potentiometer Colors: Pink, tan, cocoa, medium brown	Remote controls WatchPlot and SoundPlot FM TacTronic program selection switch On/Off switch (battery cover) Colors: Pink, tan, cocoa, medium brown	Remote controls WatchPlot and SoundPlot FM TacTronic program selection switch On/Off switch (battery cover) Colors: Pink, tan, cocoa, medium brown
Programming	Loudness Perception Profile (LPP) with Active Response Unit (ARU) Fitting formulas: Phonak Digital, Phonak Digital Ski Slope and DSL (i/o)	Loudness Perception Profile (LPP) with Active Response Unit (ARU) Fitting formulas: Phonak Digital, Phonak Digital Ski Slope and DSL (i/o)	Loudness Perception Profile (LPP) with Active Response Unit (ARU) Fitting formulas: Phonak Digital, Phonak Digital Ski Slope and DSL (i/o)

Claims	Prior Art
	<p>Lewis – Speech perception in noise: directional microphones (2004) at Page 431</p> <p>All subjects were fit with digital Phonak Claro 311 dAZ BTE hearing aids bilaterally. All earmolds had select-a-vent (SAV) venting and #13 or 3 mm horn tubing. In addition to the hearing aids, subjects were fit with Phonak Microlink ML8 FM receivers bilaterally. These FM receivers attach to the bottom of a BTE hearing aid and may be utilized in either the “FM only” mode, which attenuates the hearing aid microphone by 20 dB, or in the FM plus hearing aid mode, which allows for FM input and input of environmental sounds via the hearing aid microphones simultaneously at the same output level. The Phonak TX3 HandyMic FM transmitter served as the FM transmitter. This transmitter has three microphone options: (1) Wide Angle, which picks up sounds arriving from all directions around the transmitter microphone equally; (2) Zoom, which provides reduced amplification to signals arriving from the rear (cardioid); and (3) SuperZoom, which provides reduced amplification to signals originating from the rear and the sides (hypercardioid). The</p>

Claims	Prior Art
	<p data-bbox="583 175 1449 207">Oticon – DigiFocus Compact User Instructional Brochure at Page 4</p>  <p data-bbox="583 841 1075 873">Oticon – DigiFocus Description at p. 1</p> <p data-bbox="634 896 1312 966">DigiFocus is available both as behind-the-ear models (BTE) in colours which match the hair tone, and as in-the-ear (ITE) models which fit discreetly into the ear.</p>  <p data-bbox="714 1291 1239 1315"><i>DigiFocus Compact, Canal, Half Shell and Low Profile.</i></p>

Claims	Prior Art
	<p data-bbox="583 175 1115 207">Oticon – DigiFocus Description at Page 1</p> <p data-bbox="636 232 1316 391">Inside each DigiFocus, the digital chips have the same processing power as a fast laptop computer. This power is used to analyse the incoming sound signal, and speech is extracted through a unique system called ASA - Adaptive Speech Alignment. ASA always makes speech audible, amplifying it as much as necessary without making it uncomfortably loud or distorted. Sounds are reproduced very naturally - making speech far easier to understand than with any other hearing instrument on the market.</p> <p data-bbox="583 423 1220 456">Oticon – DigiFocus White Paper (1996) at Page 3</p> <p data-bbox="636 472 1299 680">After the A/D-conversion, the resulting digital code can be mathematically manipulated, or processed (as it is normally expressed when talking about sound). This is where the real benefits of digital technology come in. The amount and types of processing that are possible with digital technology far surpass the capabilities of analogue processing. Once the processing is complete, the digital signal is converted back to analogue form by the D/A-converter. It is now ready for transmission to the small receiver inside the hearing instrument, which makes the sound audible. □</p> <p data-bbox="583 712 1234 745">Oticon – DigiFocus White Paper (1996) at Page 11</p> <p data-bbox="636 769 1308 899">Many people find it uncomfortable to be wired to a computer. This "constricted" feeling can be avoided by using the OtiSet programming belt, which is worn on the shoulder. The programming belt provides cordless infrared communication between the PC and the hearing instrument while in the ear. Alternatively, you may use the HI-PRO interface for a wired solution.</p> <p data-bbox="636 924 1308 1079">OtiSet offers a choice of three different views of the fitting. The Audiological View provides a visual guide to instrument performance in four typical sound environments. The Technical View illustrates traditional technical information such as coupler data and insertion gain at different input levels. With the Binaural View you are able to examine settings for both right and left instruments at the same time.</p>

Claims	Prior Art
	<p data-bbox="583 175 894 207">Anderson '783 at FIG. 1</p>  <p data-bbox="995 748 1062 776">FIG. 1</p> <p data-bbox="583 846 919 878">Anderson '783 at 3:13–18</p> <p data-bbox="632 883 1969 1024">FIG. 1 is a block diagram system overview illustrating how an earpiece (comprising a microphone, RF transceiver, dual-use antenna/extractor and speaker) is worn in the ear, communicates via a primary two-way RF link with an RPU, and also communicates via an optional secondary wireless link to a telephone system.</p> <p data-bbox="583 1081 945 1114">Anderson '783 at 3:52–4:14</p> <p data-bbox="632 1118 1986 1442">One preferred embodiment of the invention herein described is shown in block diagram form in FIG. 1. In this figure, the earpiece 10 shown worn in the ear 11 uses a standard completely-in-the-canal (CIC) housing, well known in the art, although many other earpiece housing types known in the art (e.g., behind-the-ear, or BTE) can be used. The earpiece shown in FIG. 1 comprises a microphone 12, RF transceiver 13 with an antenna that doubles as an earpiece extractor 14, and speaker 15. Although the dual-use antenna/extractor arrangement shown in FIG. 1 is often desirable, a separate antenna or antennae can be used in conjunction with many different types of earpiece extractors as appropriate for a particular application. Antennae and extractors can be disguised as ear hair or jewelry (e.g., earrings), or conductive filaments can be permanently implanted in the ear cartilage to act as antennae using a process similar to</p>

Claims	Prior Art
	<p>ear piercing. The earpiece 10 communicates with a remote processor unit (RPU) 16 via a primary two-way RF link 17, although many other wireless link media (e.g., ultrasonic or infrared) may be used instead. The RPU is typically worn under clothing (e.g., carried in a pocket or purse), but may also be worn in plain sight (e.g., on a belt) if desired. The RPU 16 may be connected (via wired or wireless means 18) to optional secondary wireless link circuitry 19 that allows wireless communication between the RPU and other sources of information (e.g., the general subscriber telephone network) via a secondary wireless link. Note that the optional secondary wireless link circuitry 19 may or may not be contained within the RPU case 70 of FIG. 7.</p> <p>Anderson '783 at 4:30–33</p> <p>Although the system described here is monaural, it will be clear to those skilled in the art that a binaural (stereo) system can be constructed based on the same principles as the monaural system.</p> <p>Anderson '783 at 22:11–61</p> <p>Stand-alone earpiece operation: If an earpiece cannot communicate with an RPU for any reason, stand-alone earpiece operation as a CIC hearing aid is generally desirable. Examples include occasions where the RPU must be turned off by pushing the RPU power-off key 78 (e.g., situations where operation of wireless transmitters is prohibited such as aircraft takeoffs and classified military briefings), RPU battery or other RPU electronics failure, misplacement of an RPU, and primary wireless link 17 interference. Note that a microphone located in the earpiece 10 is preferred to a body-worn microphone, since only an earpiece microphone 12 can be used to provide the desired default CIC hearing aid function (earpiece microphones also provide superior audio localization capabilities as the microphone moves with the user's head, and provide reduced wind and clothing noise compared with body-worn microphones). Parameters that control the stand-alone earpiece operation are stored in the earpiece CIC default parameter storage 870, which is accessed by the ASIC 845 to implement the desired stand-alone functions as described earlier. CIC default parameters include, but are not limited to, settings for the nominal and maximum allowable earpiece volume. Other information contained in the default parameter storage 870 includes, but is not limited to, the communication address and L/R (left/right) selection bit. To ensure that the CIC default parameter storage information cannot be changed inadvertently, changes are only allowed when the earpiece is not in use as, for example, when the earpiece battery 860 is being charged. A battery charging detector 840, that is sensitive to the high-level magnetic induction fields commonly used for earpiece battery charging, indicates to the ASIC 845 that the default parameter storage 870 values may be changed upon receipt of an appropriate command from the RPU. Specifically, to change information in the earpiece CIC default parameter storage 870, a user first activates the earpiece battery charging system,</p>

Claims	Prior Art
	<p>e.g., by placing the earpiece in a recharger/storage tray well known in the art. The RPU keyboard 946 is then used to enter a command to the DSP 948 that indicates the user's desire to change information in the default parameter storage 870, followed by the new information to be stored. The DSP 948 causes the RPU to transmit an interrogation with a reserved address code to which all earpieces being charged within the communication range of the RPU are receptive. The information that is to ultimately be loaded into the earpiece default parameter storage 870 is then transmitted by the RPU using the auxiliary bit field of many successive interrogations, said interrogations are received by the earpiece, and the desired information is loaded into the CIC default parameter storage 870 under control of the ASIC 845.</p> <p>Anderson '783 at 22:61–23:35</p> <p>Details of normal system operation: During normal operation of the hearing aid with wireless remote processor, sounds from the ambient environment are picked up by a microphone 12 in the earpiece 10 and transmitted with other information (e.g., data representing the earpiece battery 860 voltage as obtained by an A/D 865 converter and transmitted under control of the earpiece ASIC 845 in the auxiliary bit of many successive replies via the reply data register 825) over a primary two-way wireless link 17 to the RPU 16, where the audio signals are enhanced according to the user's needs by the RPU DSP 948. Signal enhancement is achieved via any of a number of techniques well known in the art that use the capabilities of the general-purpose RPU DSP 948. The enhanced audio signals may be combined with other information (e.g., a synthesized voice, generated within the DSP 948, warning that the aforementioned earpiece battery voltage is low) and transmitted from the RPU 16 over the primary wireless link 17 to the earpiece 10, where they are converted by a speaker 15 to sounds that can only be heard by the user 11. Art auxiliary speaker can be supplied in the RPU 16 as a peripheral 950 device that allows people other than the user to hear the same sounds as those produced by the earpiece speaker 15, if desired (as, e.g., when a person other than the user wishes to use the system's optional cellular telephone capabilities). A different, optional secondary two-way wireless link (shown as part of the RPU in FIG. 9 as secondary wireless link circuitry 944, or connected to the RPU by a wired or other, e.g. infrared, wireless link 18 and shown as optional secondary wireless link circuitry 19 in FIG. 1) can be used for communication between the RPU and a cellular telephone system or other sources of information. An RPU keyboard suitable for covert data entry 72, or voice recognition capabilities in the RPU (implemented, e.g., using any of a number of techniques well known in the art that use the general-purpose RPU DSP 948), can be used to control hearing aid parameters (such as amplification level set by the RPU DSP 948) and telephone dialing functions (e.g., automatic dialing of a predefined telephone number using information stored in the RPU DSP 948) in a manner that can be made imperceptible to a casual observer if desired.</p>

Claims	Prior Art
	<p>A POSITA would have been motivated to combine the teachings of the DHA system with Anderson to meet the claim limitations with a reasonable expectation of success for the reasons described above with respect to claims [1a] and 11.</p> <p>To the extent that the DHA system, alone or in view of Anderson, does not disclose any elements of this claim limitation, it would have been obvious to a person having ordinary skill in the art to combine, with the teachings of DHA systems, those elements, as known to them or as disclosed by other prior art identified in the cover pleading.</p>
<p>14. A personal audio system as recited in claim 1, wherein the mobile communication device is or includes a mobile phone.</p>	<p>DHA systems, alone or in view of Anderson, expressly or inherently discloses “[a] personal audio system as recited in claim 1, wherein the mobile communication device is or includes a mobile phone.”</p> <p>Phonak - SmartLink SX - FM Systems Manual at p. 3</p> <p>SmartLink SX is the first communication device to incorporate</p> <ul style="list-style-type: none"> • Advanced FM features • Hearing Instrument remote control • A Bluetooth link for cell phone use <p>Phonak - SmartLink SX - FM Systems Manual at p. 10</p> <div> <div> <p>Smart links with Bluetooth</p> <p>With the inclusion of Bluetooth technology in SmartLink, Phonak now enables hearing instrument wearers to enjoy unrestricted use of all Bluetooth compatible cellular phones. One of the major complaints of hearing impaired people is now history – talking on the telephone is even easier for them than for people with normal hearing. Simply synchronizing the Bluetooth telephone and SmartLink creates the communication path. As long as the cell phone is within a 10-meter range of SmartLink, all incoming calls are automatically directed</p> </div> <div> <p>via SmartLink to the hearing instrument. The cell phone is only needed for dialing outgoing calls. Whenever the cell phone is switched on, the Bluetooth function remains active for incoming calls, regardless of whether SmartLink is switched on or off. This guarantees that telephone calls are not missed, even when the FM function is not required. When using an external audio input to listen to music, incoming telephone calls have priority when the cell phone is switched on. The audio signal is automatically decreased so that the user hears the ringing tone clearly.</p> </div> </div>

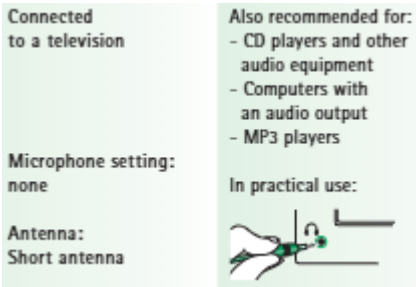
Claims	Prior Art
	<p>Phonak - SmartLink SX - FM Systems Manual at p. 11</p> <p>Smart digital synthesizer technology</p> <p>Digital frequency synthesizer technology, only available with Phonak FM communication systems, enables online frequency synchronization and scanning of the frequencies used. Once a frequency is selected on the SmartLink, it automatically sends a signal to the FM receiver that adjusts the frequency of the receiver accordingly. This automatic synchronization process occurs every time the SmartLink is switched on. In cases where more than one FM transmitter is in use, at lectures or presentations for example, the SmartLink has an exclusive scan function. This seeks out the frequency being used by the other transmitter and sets the receiver frequency accordingly.</p> <p>Phonak - SmartLink SX - FM Systems Manual at p. 14</p> 

Claims	Prior Art
	<p>Phonak - SmartLink SX User Guide at mobi</p> <ul style="list-style-type: none"> • SmartLink can communicate with a Bluetooth compatible mobile phone. Please ask your mobile phone dealer if you are uncertain whether your model is Bluetooth-capable. <p>Phonak – Microlink User Guide at Page 7</p> <p>You may even connect the system directly to your telephone, if it has an audio output socket. Alternatively, use the optional adapter available for your HandyMic. Switching off the microphone then allows you to hear the telephone signal without any background noise.</p> <p>Phonak – Annual Report (2003–2004) at p. 12</p> <p>Within the FM business segment, the Group continued to leverage existing key technology into new applications. With SmartLink, Phonak has achieved another industry first by allowing hearing instrument users to access other communications systems. SmartLink combines a remote control, a transmitter, a directional microphone system and a Bluetooth link to mobile phone all in one single device. Phonak also presented EduLink, an innovative FM device to assist the learning development of children with Auditory Processing Disorder (APD). EduLink was perceived very positively within the medical community and recent classroom tests show favourable results.</p>

Claims	Prior Art
	<p data-bbox="583 175 1226 211">Phonak SmartLink SX Product Information at p. 1</p> <div data-bbox="630 224 1129 818"> <p data-bbox="642 230 856 256">Product Information</p> <p data-bbox="642 263 982 321">SmartLink™ SX</p>  <p data-bbox="642 659 1121 818">SmartLink SX is a compact communication device for people with hearing loss who lead busy, active lives. The SmartLink SX incorporates an FM transmitter with digital signal processing, a remote control for Phonak hearing instruments, and Bluetooth for communication with mobile phones. With its digital frequency synthesizer, the SmartLink allows a highly flexible choice of operating frequency. For maximum customer benefit, we recommend using a Phonak multi-channel FM receiver (e.g. the MLxS or MLBS) that can be automatically synchronized from the SmartLink.</p> </div>

Claims	Prior Art
	<p>Phonak – Collections 2002 Product List at Page 29</p>  <p>Anderson '783 at 1:50–65</p> <p>The present invention uses an audio path that is different from the audio path used by devices in the prior art. The present invention uses a primary two-way wireless link (not a one-way link as in the prior art) between an earpiece worn in or at the user's ear and an RPU worn by or located near the user. In the present invention, audio signals from the environment are picked up by a microphone in the earpiece (not a hand-held microphone as used by wireless hearing aids in the prior art) and transmitted over the primary two-way wireless link to the RPU (instead of going to processor electronics contained in the earpiece, as in basic hearing aids known in the prior art), where the audio signals are enhanced according to the user's needs before transmission over the primary wireless link to the earpiece. Signal processing is performed in the RPU rather than the earpiece to take advantage of relaxed size and power constraints.</p> <p>Anderson '783 at 2:40–59</p> <p>Another object of this invention is to provide supplemental audio information (e.g. a verbal warning from the RPU that an earpiece battery is low) and communication services (e.g., cellular telephone and paging services) to the user via the hearing aid system. Such services are accessed in a manner that can be made imperceptible to a casual observer if desired. The RPU (which contains a digital signal processor or other</p>

Claims	Prior Art
	<p>computer) acts as one source of information, e.g., by using a synthesized voice message to provide the time of day, and the RPU can also be used to access a secondary wireless link to the general subscriber telephone network or voice paging services. User control of hearing aid parameters and requests for information are accomplished using pushbuttons located on the RPU (including pushbuttons suitable for data entry in a covert manner) or voice recognition. Many sensors and peripheral devices can reside in or be attached to the RPU by wired or wireless means, and can provide a variety of information for different applications (e.g., heart pulse rate) as audio in the user's earpiece.</p> <p>Anderson '783 at 24:47–66</p> <p>Morse code keys 956 or voice recognition capabilities implemented using the RPU DSP 948) in a manner that can be made imperceptible to a casual observer if desired. In this respect, the system acts as a "virtual handset" communication aid that replaces a conventional telephone handset in many applications (although users may additionally communicate with a conventional handset or speakerphone in the usual fashion). To accomplish the virtual handset function, the RPU 16 is equipped with optional secondary two-way wireless link circuitry 19 to communicate with a telephone or other communication network. The secondary wireless link is separate from, and does not interfere with, the primary wireless link 17 between the earpiece 10 and RPU 16. In applications where relatively high-power transmissions take place on the secondary link and the secondary link transmitter is integrated into the RPU, the RPU can be placed on a user's wrist or ankle away from vital organs to reduce any potential health risk relative to conventional wireless handsets.</p> <p>A POSITA would have been motivated to combine the teachings of the DHA system with Anderson to meet the claim limitations with a reasonable expectation of success for the reasons described above with respect to claim 1b.</p> <p>To the extent that the DHA system, alone or in view of Anderson, does not disclose any elements of this claim limitation, it would have been obvious to a person having ordinary skill in the art to combine, with the teachings of DHA systems, those elements, as known to them or as disclosed by other prior art identified in the cover pleading.</p>
<p>15. A personal audio system as recited in claim 1, wherein the customized audio output includes music.</p>	<p>DHA systems, alone or in view of Anderson, expressly or inherently discloses “[a] personal audio system as recited in claim 1, wherein the customized audio output includes music.”</p>

Claims	Prior Art
	<p>Phonak - SmartLink SX User Guide at p. 27</p> <ul style="list-style-type: none"> • The SmartLink can be directly connected to audio equipment like televisions, CD and MP3 players, and computers with audio outputs. <p>Phonak - SmartLink SX User Guide at p. 47</p>  <p>The screenshot shows a user interface for the SmartLink SX device. It is divided into two columns. The left column has a light green background and contains the following text: 'Connected to a television', 'Microphone setting: none', and 'Antenna: Short antenna'. The right column has a light green background and contains the following text: 'Also recommended for:' followed by a bulleted list: '- CD players and other audio equipment', '- Computers with an audio output', and '- MP3 players'. Below this list, it says 'In practical use:' followed by a small diagram showing a hand holding a device with a green antenna pointing towards a television.</p> <p>Anderson '783 at 1:50–65</p> <p>The present invention uses an audio path that is different from the audio path used by devices in the prior art. The present invention uses a primary two-way wireless link (not a one-way link as in the prior art) between an earpiece worn in or at the user's ear and an RPU worn by or located near the user. In the present invention, audio signals from the environment are picked up by a microphone in the earpiece (not a hand-held microphone as used by wireless hearing aids in the prior art) and transmitted over the primary two-way wireless link to the RPU (instead of going to processor electronics contained in the earpiece, as in basic hearing aids known in the prior art), where the audio signals are enhanced according to the user's needs before transmission over the primary wireless link to the earpiece. Signal processing is performed in the RPU rather than the earpiece to take advantage of relaxed size and power constraints.</p> <p>Anderson '783 at 2:40–59</p> <p>Another object of this invention is to provide supplemental audio information (e.g. a verbal warning from the RPU that an earpiece battery is low) and communication services (e.g., cellular telephone and paging services) to the user via the hearing aid system. Such services are accessed in a manner that can be made imperceptible to a casual observer if desired. The RPU (which contains a digital signal processor or other computer) acts as one source of information, e.g., by using a synthesized voice message to provide the time of day, and the RPU can also be used to access a secondary wireless link to the general subscriber telephone network or voice paging services. User control of hearing aid parameters and requests for information are accomplished using pushbuttons located on the RPU (including pushbuttons suitable for data entry in a covert manner) or voice recognition. Many sensors and peripheral devices can reside in or</p>

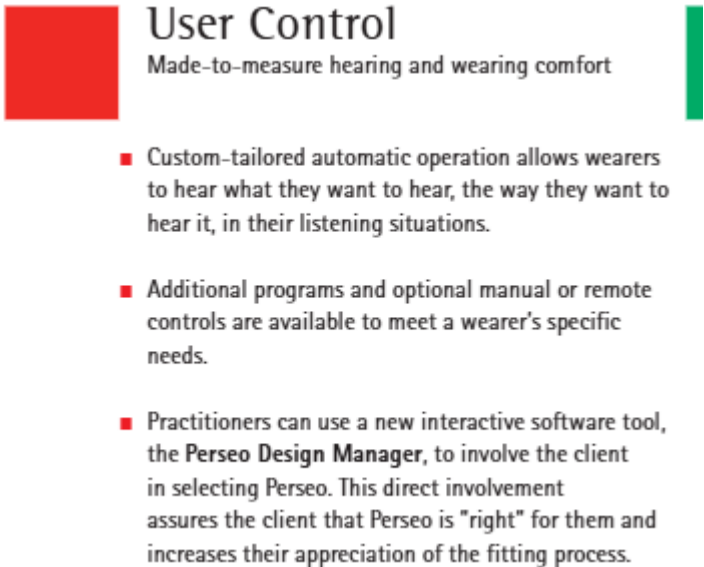
Claims	Prior Art
	<p>be attached to the RPU by wired or wireless means, and can provide a variety of information for different applications (e.g., heart pulse rate) as audio in the user's earpiece.</p> <p>Anderson '783 at 24:15–26 The invention provides useful unsolicited information (i.e., information automatically provided to the user by the system, not as a result of a user request for information) to the user of the system in a manner that can be made imperceptible to a casual observer if desired. Information about the status of the hearing aid system itself, such as a verbal warning that an earpiece or RPU battery voltage is low, is provided to the user as detailed earlier. Automatic measurements of heart pulse rate, body temperature and walking distance traveled can be taken by sensors 950 in the RPU 16 or earpiece 10 and stored in the RPU DSP 948.</p> <p>Anderson '783 at 24:26–32 Sensor measurements can be reported to the user by a voice synthesized in the RPU DSP 948 using any of a variety of synthesis techniques well known in the art, including playback of pre-recorded digitized speech. The voice signal produced within the RPU DSP 948 is then added to the processed audio signals already being sent from the RPU to the earpiece.</p> <p>Anderson '783 at 24:45–47 Solicited information: Solicited information is provided by the system (e.g., via the earpiece speaker 15) in response to user requests (e.g., via an RPU keyboard 946,</p> <p>Anderson '783 at 24:47–66 Morse code keys 956 or voice recognition capabilities implemented using the RPU DSP 948) in a manner that can be made imperceptible to a casual observer if desired. In this respect, the system acts as a "virtual handset" communication aid that replaces a conventional telephone handset in many applications (although users may additionally communicate with a conventional handset or speakerphone in the usual fashion). To accomplish the virtual handset function, the RPU 16 is equipped with optional secondary two-way wireless link circuitry 19 to communicate with a telephone or other communication network. The secondary wireless link is separate from, and does not interfere with, the primary wireless link 17 between the earpiece 10 and RPU 16. In applications where relatively high-power transmissions take place on the secondary link and the secondary link transmitter is integrated into the RPU, the RPU can be placed on a user's wrist or ankle away from vital organs to reduce any potential health risk relative to</p>

Claims	Prior Art
	<p>conventional wireless handsets.</p> <p>Anderson '783 at 25:50–26:5</p> <p>Sensors and peripheral devices: Many sensors and peripheral devices 950 can reside in or be attached to the RPU by wired or wireless means, and controlled by the user (via inputs to the RPU DSP 948 through the keyboard 946, Morse code keys 956 or voice recognition) to provide both solicited and unsolicited information. Music can be provided from pre-stored information including tape, compact disc, over-the-telephone services via the secondary wireless link, or a commercial broadcast radio receiver. User location information can be provided to the RPU from special-purpose transmitters at specific locations in nursing home applications. Self-contained user-carried inertial navigation packages (with solid-state gyroscopes, accelerometers and a compass) or a Global Positioning System receiver can also be used. Sensors for a foot soldier RPU application may include many types of detectors such as RF detectors (e.g., radar detection capability), nuclear/biological/chemical weapon detectors, infrared or laser detectors and metal detectors. Appropriate verbal warnings or distinctive alerting tones generated by the RPU DSP 948 are sent from the RPU 16 to the earpiece 10 and ultimately received as audio signals by the user.</p> <p>A POSITA would have been motivated to combine the teachings of the DHA system with Anderson to meet the claim limitations with a reasonable expectation of success for the reasons described above with respect to claim [1b].</p> <p>To the extent that the DHA system, alone or in view of Anderson, does not disclose any elements of this claim limitation, it would have been obvious to a person having ordinary skill in the art to combine, with the teachings of DHA systems, those elements, as known to them or as disclosed by other prior art identified in the cover pleading.</p>
<p>16. A personal audio system as recited in claim 1, wherein the customized audio output includes a voice conversation.</p>	<p>DHA systems, alone or in view of Anderson, expressly or inherently discloses “[a] personal audio system as recited in claim 1, wherein the customized audio output includes a voice conversation.”</p> <p>DHA Systems are designed to aid a user’s hearing, and a POSITA would have recognized that voice conversations were one of the expected inputs to a hearing aid, as people engage in such conversations on a daily basis.</p> <p>Anderson '783 at 5:3–21</p> <p>During normal operation, speech signals from a nearby talker (and other signals in the ambient audio environment) are picked up by the earpiece wireless microphone 20 and transmitted to the RPU FM receiver 24. The RPU FM receiver 24 output level may be adjusted using the RPU FM receiver 24 volume</p>

Claims	Prior Art
	<p>control. The resulting electrical waveform representing signals from the nearby talker travels through the DPDT switch 25, which is set to the lower position as shown in FIG. 2, to the signal enhancer 26. The signal enhancer 26 may be, for example, a voice changer device that varies the pitch of a received speech signal according to the setting of pushbutton controls located on the RPU signal enhancer 26. The speech signal's pitch can then be raised or lowered as desired to help compensate for a user's hearing loss relative to a particular talker's voice characteristics. The modified speech signal travels from the RPU signal enhancer 26 to the RPU FM transmitter 27, and finally to the earpiece headset 21 receiver where the signal is converted to acoustic waves heard by the user.</p> <p>A POSITA would have been motivated to combine the teachings of the DHA system with Anderson, to meet the claim limitations with a reasonable expectation of success for the reasons described above with respect to claim 1b.</p> <p>To the extent that the DHA system, alone or in view of Anderson, does not disclose any elements of this claim limitation, it would have been obvious to a person having ordinary skill in the art to combine, with the teachings of DHA systems, those elements, as known to them or as disclosed by other prior art identified in the cover pleading.</p>
<p>[17a] A personal audio system as recited in claim 1, wherein at least one of the earbuds in the pair of earbuds comprises a plurality of microphones configured for audio pickup,</p>	<p>DHA systems, alone or in view of Anderson, expressly or inherently discloses “[a] personal audio system as recited in claim 1, wherein at least one of the earbuds in the pair of earbuds comprises a plurality of microphones configured for audio pickup.”</p> <p><i>See above at claim 6.</i></p> <p>A POSITA would have been motivated to combine the teachings of the DHA system with Anderson to meet the claim limitations with a reasonable expectation of success for the reasons described above with respect to claim [1a].</p> <p>To the extent that the DHA system, alone or in view of Anderson, does not disclose any elements of this claim limitation, it would have been obvious to a person having ordinary skill in the art to combine, with the teachings of DHA systems, those elements, as known to them or as disclosed by other prior art identified in the cover pleading.</p>
<p>[17b] wherein at least one of the microphones is configured to pick up an audio command, and</p>	<p>DHA systems, alone or in view of Anderson, expressly or inherently discloses “wherein at least one of the microphones is configured to pick up an audio command.”</p>

Claims	Prior Art
	<p>Anderson '783 at 25:29–30 The verbal commands are received as audio signals by the earpiece 10 and transmitted to the RPU 16.</p> <p>A POSITA would have been motivated to combine the teachings of the DHA system with Anderson to meet the claim limitations with a reasonable expectation of success for the reasons described above with respect to claim [1a].</p> <p>To the extent that the DHA system, alone or in view of Anderson, does not disclose any elements of this claim limitation, it would have been obvious to a person having ordinary skill in the art to combine, with the teachings of DHA systems, those elements, as known to them or as disclosed by other prior art identified in the cover pleading.</p>
<p>[17c] wherein the mobile communication device is configured to induce an operation or action based on the audio command.</p>	<p>DHA systems, alone or in view of Anderson or Seshadri, expressly or inherently discloses “wherein the mobile communication device is configured to induce an operation or action based on the audio command.”</p> <p>Anderson '783 at 25:15–27 To obtain solicited information, user requests are processed in a variety of means well known in the art, including inputs via one or more pushbuttons on the RPU and voice recognition computer programs implemented in the RPU electronics. Voice recognition can be used almost exclusively, if desired, thereby eliminating nearly all RPU pushbutton controls. Voice commands recognized by the RPU can cause the RPU to initiate actions such as changing hearing aid parameter settings (e.g., volume), or activate appliances and environmental controls (e.g., aids for the physically handicapped) in a distributed wireless network via the secondary wireless link.</p> <p>Anderson '783 at 25:27–29 RPU circuitry may also enable a user to voice-dial a desired telephone number by speaking certain pre-defined commands.</p> <p>Seshadri '839 at [0042] Input module 128, which may be a keypad, touch screen, voice recognition circuit, or other like user interfaces, receives user commands and produces digital command messages 124 there from. Such digital command messages 124 includes, but are not limited to, packet size, synchronization information, frequency hopping initiation information, timeslot allocation information, link establishment information, piconet address information, fast-forward, play, pause, volume adjust, record, stop and rewind.</p>

Claims	Prior Art
	<p>Seshadri '839 at [0047]</p> <p>Detachable ear-piece 156 and microphone 154 may have on-chip operations to support call conferencing, call waiting, flash, and other features associated with telephones. These functions may be accessed and reviewed by a user interface 158 and display 160 within host device 152 or a user interface and display located on either detachable ear-piece 156 or microphone 154. The user interface and display, located on either the host device or detachable ear-piece 156 and microphone 154 may have a display and button(s) that may be used to program device, perform directory functions including selecting number to call, view caller ID, initiate call waiting, or initiate call conferencing. Additionally, circuitry within the ear-piece 156 and microphone 154 may enable voice activated dialing. The actual voice recognition could be performed within ear-piece 156, microphone 154, or host device 152. Thus, ear-piece 156 and microphone 154 may act to initiate calls and receive calls.</p> <p>A POSITA would have been motivated to combine the teachings of the DHA system with Anderson or Seshadri to meet the claim limitations with a reasonable expectation of success for the reasons described above with respect to claim [1a].</p> <p>To the extent that the DHA system, alone or in view of Anderson, does not disclose any elements of this claim limitation, it would have been obvious to a person having ordinary skill in the art to combine, with the teachings of DHA systems, those elements, as known to them or as disclosed by other prior art identified in the cover pleading.</p>
<p>18. A personal audio system as recited in claim 1, wherein the customized audio output signals are produced based on amplifying or attenuating different frequency bands differently according to at least the stored hearing profile.</p>	<p>DHA systems, alone or in view of Anderson, expressly or inherently discloses “[a] personal audio system as recited in claim 1, wherein the customized audio output signals are produced based on amplifying or attenuating different frequency bands differently according to at least the stored hearing profile.”</p> <p><i>See</i> above at claim [1d].</p> <p>A POSITA would have been motivated to combine the teachings of the DHA system with Anderson to meet the claim limitations with a reasonable expectation of success for the reasons described above with respect to claim [1b].</p> <p>To the extent that the DHA system, alone or in view of Anderson, does not disclose any elements of this claim limitation, it would have been obvious to a person having ordinary skill in the art to combine, with the teachings of DHA systems, those elements, as known to them or as disclosed by other prior art identified in the cover pleading.</p>

Claims	Prior Art
<p>19. A personal audio system as recited in claim 1, wherein the hearing profile associated with the user is self-determined by the user.</p>	<p>DHA systems, alone or in view of Anderson, expressly or inherently discloses “[a] personal audio system as recited in claim 1, wherein the hearing profile associated with the user is self-determined by the user.”</p> <p>Phonak – Perseo Brochure at Page 7</p>  <p>The brochure titled "User Control" features a red square on the left and a green vertical bar on the right. The text reads: "User Control Made-to-measure hearing and wearing comfort". It lists three bullet points: "Custom-tailored automatic operation allows wearers to hear what they want to hear, the way they want to hear it, in their listening situations.", "Additional programs and optional manual or remote controls are available to meet a wearer's specific needs.", and "Practitioners can use a new interactive software tool, the Perseo Design Manager, to involve the client in selecting Perseo. This direct involvement assures the client that Perseo is 'right' for them and increases their appreciation of the fitting process."</p> <p>Phonak – Fitting Guide for MLxi and MLxi Baha at Page 4</p>  <p>The guide section is titled "4.2 Advanced fitting of MLxi using FM SuccessWare". It states: "The extended mode of FM SuccessWare enables you to fine-tune some parameters of the MLxi according to the needs of your customers." Below this is a screenshot of the FM SuccessWare software interface. The interface shows a sidebar with "Extended Fitting" selected, and a main area with "Personal preferences" and "Associated hearing instrument" settings. A "Basic Type" dropdown is set to "200 Hz normal". A "Gain" slider is positioned between -0.40 and 0.40. Other options include "Binaural", "Unilateral", "Bilateral", and "Unilateral".</p>

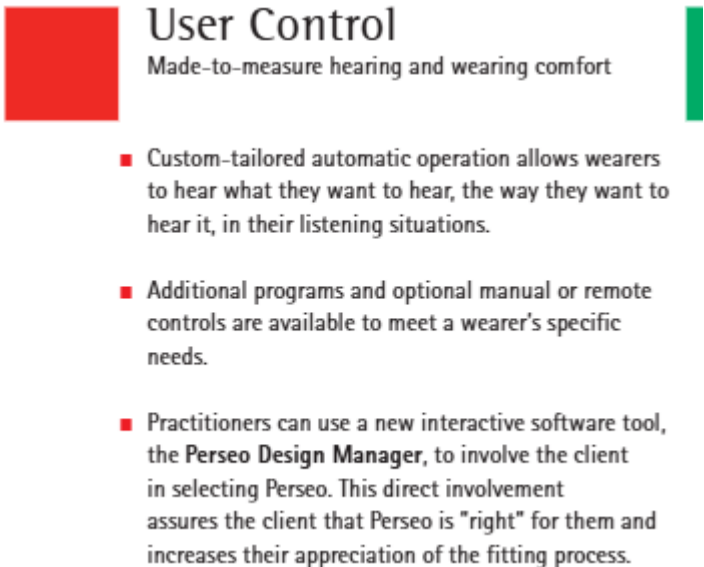
Claims	Prior Art
	<p>Phonak – Fitting Guide for MLxi and MLxi Baha at Page 8</p> <div data-bbox="636 232 911 264">7. Programming options</div> <div data-bbox="636 318 856 342">7.1 Programming interface</div> <div data-bbox="690 363 1110 607">  <p>The diagram shows a white rectangular device labeled 'FM Programming Interface' and 'PHONAK'. A black cable is plugged into the bottom of the device. A hearing aid is shown with a small black cable connected to it, which is then plugged into the top of the FM Programming Interface device.</p> </div> <div data-bbox="636 631 1110 776"> <p>All receivers can be programmed using the FM programming interface:</p> <ol style="list-style-type: none"> 1. Make sure all equipment is charged or a fresh battery is inserted in the hearing instrument. Turn all equipment on. 2. Connect the universal receiver to the hearing instrument and insert it into the programming interface. 3. Always select the option "via the programming interface" when FM SuccessWare requests communication with the receiver. </div> <div data-bbox="1129 237 1274 261">7.2 Using inspiro</div> <div data-bbox="1150 272 1604 630">  <p>The diagram shows a small white device labeled 'inspiro' connected to a hearing aid by a black cable. Another black cable connects the inspiro device to a laptop.</p> </div> <div data-bbox="1129 631 1604 776"> <p>MLxi, MLxi Baha can be programmed via inspiro:</p> <ol style="list-style-type: none"> 1. Make sure a fresh battery is inserted in the hearing instrument. Turn all equipment on. 2. Connect the universal receiver to the hearing instrument, the CI speech processor, or MicroLink CI S, and place in front of the inspiro. 3. Always select the option "via inspiro" when the FM SuccessWare requests communication with the receiver. </div> <p>Phonak – Media Release re Product Offerings in F2003/2004 at Page 2</p> <p>The key breakthrough technology in Perseo is PersonalLogic, the first ever automatic control system that is fitter programmable according to the tastes and preferences of each individual wearer. Automatic performance has now a new level of intelligence beyond merely adapting to acoustic changes in the environment. Perseo is the first hearing system with customized automatic operation, as opposed to operation based on factory settings. It allows hearing professionals to provide clients with a hearing system "custom-tailored" to their personal hearing requirements and tastes.</p>

Claims	Prior Art
	<p>Phonak – Collections 2002 Product List at Page 33</p> <div data-bbox="636 256 802 337"> <p>PFG – new perspectives in fitting success</p> </div> <div data-bbox="636 386 884 570">  </div> <div data-bbox="636 591 884 774">  </div> <div data-bbox="945 228 1297 818"> <p>QuickGuide The versatile, easy-to-use online tour through the PFG fitting process with access to all optional fitting tools.</p> <p>Instrument selection wizard This online assistant recommends the product that best matches the customer's needs.</p> <p>Initial SoundCheck This is an ideal tool for fast and simple fine tuning of sound quality. Three controls provide convenient gain and balance settings in the high- and low-frequency ranges.</p> <p>Automatic fine tuning Based on the customer's hearing perception, the software automatically computes the most suitable modifications.</p> <p>Manual fine tuning A high degree of fitting precision can be attained with this broad spectrum of parameters in an optimised display format.</p> <p>Phonak Sound The degree of fitting success can be dramatically improved by using these real-life sound clips. All sound examples can be saved on a computer as MP3 files; they are also available on a sound CD.</p> <p>Proven PFG benefits</p> <ul style="list-style-type: none"> • One program for all Phonak hearing computers • NOAH compatible or standalone version • Detailed help system with in-depth scientific and technical information • PFG training demo </div> <p>Anderson '783 at 27:4–24</p> <p>Hearing test: The present invention allows hearing tests to be performed via the earpiece normally worn by the user, without the need for additional expensive test equipment. All signal generation capabilities necessary to perform hearing tests are available in the hearing aid system comprising an earpiece 10 and RPU 16. Specifically, the RPU DSP 948 can be used to generate tones of varying frequencies and amplitudes, as well as other signals known to be useful for hearing test purposes (e.g., nonsense syllables) that are converted to acoustic waves by the earpiece. Such tones are used to perform a hearing test that determines appropriate gain vs. frequency parameters for a program stored in the RPU DSP 948 that performs signal enhancement to compensate for the user's hearing loss. Note that the hearing test program that controls the RPU DSP 948 during the hearing test can be temporarily stored in the RPU DSP 948 for the duration of the test, then deleted upon completion of the test to allow re-use of RPU DSP 948 memory resources during normal operation. The hearing test program may be loaded into the RPU DSP 948 through the secondary wireless link 944 or a wired peripheral link 950.</p>

Claims	Prior Art
	<p>A POSITA would have been motivated to combine the teachings of the DHA system with Anderson to meet the claim limitations with a reasonable expectation of success for the reasons described above with respect to claim [1b].</p> <p>To the extent that the DHA system, alone or in view of Anderson, does not disclose any elements of this claim limitation, it would have been obvious to a person having ordinary skill in the art to combine, with the teachings of DHA systems, those elements, as known to them or as disclosed by other prior art identified in the cover pleading.</p>
<p>20[pre]. A personal audio system as recited in claim 1,</p>	<p>DHA systems, alone or in view of Anderson, expressly or inherently discloses “[a] personal audio system as recited in claim 1.”</p> <p><i>See above at claim 1[pre].</i></p> <p>A POSITA would have been motivated to combine the teachings of the DHA system with Anderson to meet the claim limitations with a reasonable expectation of success for the reasons described above with respect to claim [1a] and [1b].</p> <p>To the extent that the DHA system, alone or in view of Anderson, does not disclose any elements of this claim limitation, it would have been obvious to a person having ordinary skill in the art to combine, with the teachings of DHA systems, those elements, as known to them or as disclosed by other prior art identified in the cover pleading.</p>
<p>[20a] wherein the mobile communication device includes computing capabilities and a storage device, the storage device storing at least a software program executed by the mobile communication device to enable the user to establish the hearing profile associated with the user, and</p>	<p>DHA systems, alone or in view of Anderson and Eaton, expressly or inherently discloses “wherein the mobile communication device includes computing capabilities and a storage device, the storage device storing at least a software program executed by the mobile communication device to enable the user to establish the hearing profile associated with the user.”</p> <p>Phonak - SmartLink SX - FM Systems Manual at p. 3</p> <p>SmartLink SX is the first communication device to incorporate</p> <ul style="list-style-type: none"> • Advanced FM features • Hearing Instrument remote control • A Bluetooth link for cell phone use

Claims	Prior Art
	<p data-bbox="583 175 1287 207">Phonak - SmartLink SX - FM Systems Manual at p. 10</p> <div data-bbox="636 220 1005 253"> <p>Smart links with Bluetooth</p> </div> <div data-bbox="636 259 1092 537"> <p>With the inclusion of Bluetooth technology in SmartLink, Phonak now enables hearing instrument wearers to enjoy unrestricted use of all Bluetooth compatible cellular phones. One of the major complaints of hearing impaired people is now history – talking on the telephone is even easier for them than for people with normal hearing. Simply synchronizing the Bluetooth telephone and SmartLink creates the communication path. As long as the cell phone is within a 10-meter range of SmartLink, all incoming calls are automatically directed</p> </div> <div data-bbox="1144 232 1600 537"> <p>via SmartLink to the hearing instrument. The cell phone is only needed for dialing outgoing calls. Whenever the cell phone is switched on, the Bluetooth function remains active for incoming calls, regardless of whether SmartLink is switched on or off. This guarantees that telephone calls are not missed, even when the FM function is not required. When using an external audio input to listen to music, incoming telephone calls have priority when the cell phone is switched on. The audio signal is automatically decreased so that the user hears the ringing tone clearly.</p> </div> <p data-bbox="583 565 1283 597">Phonak - SmartLink SX - FM Systems Manual at p. 11</p> <div data-bbox="636 613 978 695"> <p>Smart digital synthesizer technology</p> </div> <div data-bbox="636 698 1092 1057"> <p>Digital frequency synthesizer technology, only available with Phonak FM communication systems, enables online frequency synchronization and scanning of the frequencies used. Once a frequency is selected on the SmartLink, it automatically sends a signal to the FM receiver that adjusts the frequency of the receiver accordingly. This automatic synchronization process occurs every time the SmartLink is switched on. In cases where more than one FM transmitter is in use, at lectures or presentations for example, the SmartLink has an exclusive scan function. This seeks out the frequency being used by the other transmitter and sets the receiver frequency accordingly.</p> </div>

Claims	Prior Art
	<p data-bbox="583 175 1289 207">Phonak - SmartLink SX - FM Systems Manual at p. 14</p> <div data-bbox="630 224 1289 824">  <p>The image shows a blue Phonak SmartLink SX device at the top, which is a small electronic device with a screen and buttons. Below it is the cover of the SmartLink SX manual, which is dark blue with white text. The title 'SmartLink™ SX' is prominently displayed, followed by the tagline 'The gateway to communication'. A list of features is provided below the tagline.</p> <ul style="list-style-type: none"> • The only FM system with digital synthesizer • Cellular phone access via Bluetooth • Integrated remote control for hearing instruments • Ergonomic design with visual controls • State-of-the-art energy management </div> <p data-bbox="583 841 1150 873">Phonak - SmartLink SX User Guide at mobi</p> <ul style="list-style-type: none"> • SmartLink can communicate with a Bluetooth compatible mobile phone. Please ask your mobile phone dealer if you are uncertain whether your model is Bluetooth-capable.

Claims	Prior Art
	<p data-bbox="583 175 1050 207">Phonak – Perseo Brochure at Page 7</p> <div data-bbox="636 228 1335 792">  <p data-bbox="779 228 1245 305">User Control Made-to-measure hearing and wearing comfort</p> <ul data-bbox="751 370 1304 792" style="list-style-type: none"> <li data-bbox="751 370 1304 467">■ Custom-tailored automatic operation allows wearers to hear what they want to hear, the way they want to hear it, in their listening situations. <li data-bbox="751 500 1304 597">■ Additional programs and optional manual or remote controls are available to meet a wearer's specific needs. <li data-bbox="751 630 1304 792">■ Practitioners can use a new interactive software tool, the Perseo Design Manager, to involve the client in selecting Perseo. This direct involvement assures the client that Perseo is "right" for them and increases their appreciation of the fitting process. </div> <p data-bbox="583 816 1627 849">Valente - Effect of type of noise and loudspeaker array on the performance at 402</p> <p data-bbox="636 870 1020 1182">The Phonak fitting software (Version 8.3) was used to program the hearing aids. First, hearing thresholds were entered into the NOAH software. The manufacturer's prescribed "First-Fit" was calculated from these hearing thresholds. The default settings in each of the software "personal system" menus (provides adjustments to the "First-Fit" for prior experience; occlusion effect and feedback) were utilized. If the patient reported occlusion, then this feature was activated. If feedback was noticed, the feedback manager was accessed. All hearing</p> <p data-bbox="583 1206 1119 1239">Phonak – Microlink User Guide at Page 7</p> <p data-bbox="636 1260 1314 1365">You may even connect the system directly to your telephone, if it has an audio output socket. Alternatively, use the optional adapter available for your HandyMic. Switching off the microphone then allows you to hear the telephone signal without any background noise.</p>

Claims	Prior Art
	<p>Phonak – Annual Report (2003–2004) at p. 12</p> <p>Within the FM business segment, the Group continued to leverage existing key technology into new applications. With SmartLink, Phonak has achieved another industry first by allowing hearing instrument users to access other communications systems. SmartLink combines a remote control, a transmitter, a directional microphone system and a Bluetooth link to mobile phone all in one single device. Phonak also presented EduLink, an innovative FM device to assist the learning development of children with Auditory Processing Disorder (APD). EduLink was perceived very positively within the medical community and recent classroom tests show favourable results.</p> <p>Phonak – Fitting Guide for MLxi and MLxi Baha at Page 4</p> <p>4.2 Advanced fitting of MLxi using FM SuccessWare</p> <p>The extended mode of FM SuccessWare enables you to fine-tune some parameters of the MLxi according to the needs of your customers.</p> 

Claims	Prior Art
	<p data-bbox="583 175 1339 207">Phonak – Fitting Guide for MLxi and MLxi Baha at Page 8</p> <div data-bbox="634 230 1606 776"> <div data-bbox="634 230 1113 263">7. Programming options</div> <div data-bbox="634 321 856 344">7.1 Programming interface</div> <div data-bbox="688 370 1113 613">  <p>The diagram shows a white rectangular device labeled 'FM Programming Interface' and 'PHONAK'. A black cable is plugged into the bottom. A hearing aid is shown with a small antenna connected to the top of the device.</p> </div> <div data-bbox="634 636 1113 776"> <p>All receivers can be programmed using the FM programming interface:</p> <ol style="list-style-type: none"> 1. Make sure all equipment is charged or a fresh battery is inserted in the hearing instrument. Turn all equipment on. 2. Connect the universal receiver to the hearing instrument and insert it into the programming interface. 3. Always select the option "via the programming interface" when FM SuccessWare requests communication with the receiver. </div> <div data-bbox="1129 230 1276 263">7.2 Using inspiro</div> <div data-bbox="1150 279 1606 630">  <p>The diagram shows a small white device labeled 'inspiro' connected to a hearing aid by a thin wire. A black cable also connects the inspiro device to a laptop.</p> </div> <div data-bbox="1129 636 1606 776"> <p>MLxi, MLxi Baha can be programmed via inspiro:</p> <ol style="list-style-type: none"> 1. Make sure a fresh battery is inserted in the hearing instrument. Turn all equipment on. 2. Connect the universal receiver to the hearing instrument, the CI speech processor, or MicroLink CI S, and place in front of the inspiro. 3. Always select the option "via inspiro" when the FM SuccessWare requests communication with the receiver. </div> </div> <p data-bbox="583 808 1222 841">Phonak SmartLink SX Product Information at p. 1</p> <div data-bbox="634 857 1129 1448"> <p data-bbox="646 857 856 880">Product Information</p> <h1 data-bbox="646 896 982 945">SmartLink™ SX</h1> <div data-bbox="634 1036 1129 1263">  <p>The image shows a small, white, rectangular device with a green light on the front. It has a small antenna on the top. Below the device is a green bar with the CE mark and a small '1' in a circle.</p> </div> <p data-bbox="646 1286 1129 1448">SmartLink SX is a compact communication device for people with hearing loss who lead busy, active lives. The SmartLink SX incorporates an FM transmitter with digital signal processing, a remote control for Phonak hearing instruments, and Bluetooth for communication with mobile phones. With its digital frequency synthesizer, the SmartLink allows a highly flexible choice of operating frequency. For maximum customer benefit, we recommend using a Phonak multi-channel FM receiver (e.g. the MLxS or MLBS) that can be automatically synchronized from the SmartLink.</p> </div>

Claims	Prior Art
	<p>Phonak – Media Release re Product Offerings in F2003/2004 at Page 2</p> <p>The key breakthrough technology in Perseo is PersonalLogic, the first ever automatic control system that is fitter programmable according to the tastes and preferences of each individual wearer. Automatic performance has now a new level of intelligence beyond merely adapting to acoustic changes in the environment. Perseo is the first hearing system with customized automatic operation, as opposed to operation based on factory settings. It allows hearing professionals to provide clients with a hearing system "custom-tailored" to their personal hearing requirements and tastes.</p> <p>Phonak – Collections 2002 Product List at Page 29</p> 

Claims	Prior Art
	<p>Phonak – Collections 2002 Product List at Page 33</p> <div data-bbox="638 256 802 337"> <p>PFG – new perspectives in fitting success</p> </div> <div data-bbox="638 386 884 570">  </div> <div data-bbox="638 591 884 774">  </div> <div data-bbox="945 228 1297 818"> <p>QuickGuide The versatile, easy-to-use online tour through the PFG fitting process with access to all optional fitting tools.</p> <p>Instrument selection wizard This online assistant recommends the product that best matches the customer's needs.</p> <p>Initial SoundCheck This is an ideal tool for fast and simple fine tuning of sound quality. Three controls provide convenient gain and balance settings in the high- and low-frequency ranges.</p> <p>Automatic fine tuning Based on the customer's hearing perception, the software automatically computes the most suitable modifications.</p> <p>Manual fine tuning A high degree of fitting precision can be attained with this broad spectrum of parameters in an optimised display format.</p> <p>Phonak Sound The degree of fitting success can be dramatically improved by using these real-life sound clips. All sound examples can be saved on a computer as MP3 files; they are also available on a sound CD.</p> <p>Proven PFG benefits</p> <ul style="list-style-type: none"> • One program for all Phonak hearing computers • NOAH compatible or standalone version • Detailed help system with in-depth scientific and technical information • PFG training demo </div> <p>Anderson '783 at 1:50–65</p> <p>The present invention uses an audio path that is different from the audio path used by devices in the prior art. The present invention uses a primary two-way wireless link (not a one-way link as in the prior art) between an earpiece worn in or at the user's ear and an RPU worn by or located near the user. In the present invention, audio signals from the environment are picked up by a microphone in the earpiece (not a hand-held microphone as used by wireless hearing aids in the prior art) and transmitted over the primary two-way wireless link to the RPU (instead of going to processor electronics contained in the earpiece, as in basic hearing aids known in the prior art), where the audio signals are enhanced according to the user's needs before transmission over the primary wireless link to the earpiece. Signal processing is performed in the RPU rather than the earpiece to take advantage of relaxed size and power constraints.</p> <p>Anderson '783 at 2:40–59</p> <p>Another object of this invention is to provide supplemental audio information (e.g. a verbal warning from the RPU that an earpiece battery is low) and communication services (e.g., cellular telephone and paging services) to the user via the hearing aid system. Such services are accessed in a manner that can be made imperceptible to a casual observer if desired. The RPU (which contains a digital signal processor or other</p>

Claims	Prior Art
	<p>computer) acts as one source of information, e.g., by using a synthesized voice message to provide the time of day, and the RPU can also be used to access a secondary wireless link to the general subscriber telephone network or voice paging services. User control of hearing aid parameters and requests for information are accomplished using pushbuttons located on the RPU (including pushbuttons suitable for data entry in a covert manner) or voice recognition. Many sensors and peripheral devices can reside in or be attached to the RPU by wired or wireless means, and can provide a variety of information for different applications (e.g., heart pulse rate) as audio in the user's earpiece.</p> <p>Anderson '783 at 24:47–66</p> <p>Morse code keys 956 or voice recognition capabilities implemented using the RPU DSP 948) in a manner that can be made imperceptible to a casual observer if desired. In this respect, the system acts as a "virtual handset" communication aid that replaces a conventional telephone handset in many applications (although users may additionally communicate with a conventional handset or speakerphone in the usual fashion). To accomplish the virtual handset function, the RPU 16 is equipped with optional secondary two-way wireless link circuitry 19 to communicate with a telephone or other communication network. The secondary wireless link is separate from, and does not interfere with, the primary wireless link 17 between the earpiece 10 and RPU 16. In applications where relatively high-power transmissions take place on the secondary link and the secondary link transmitter is integrated into the RPU, the RPU can be placed on a user's wrist or ankle away from vital organs to reduce any potential health risk relative to conventional wireless handsets.</p> <p>Anderson '783 at 27:4–24</p> <p>Hearing test: The present invention allows hearing tests to be performed via the earpiece normally worn by the user, without the need for additional expensive test equipment. All signal generation capabilities necessary to perform hearing tests are available in the hearing aid system comprising an earpiece 10 and RPU 16. Specifically, the RPU DSP 948 can be used to generate tones of varying frequencies and amplitudes, as well as other signals known to be useful for hearing test purposes (e.g., nonsense syllables) that are converted to acoustic waves by the earpiece. Such tones are used to perform a hearing test that determines appropriate gain vs. frequency parameters for a program stored in the RPU DSP 948 that performs signal enhancement to compensate for the user's hearing loss. Note that the hearing test program that controls the RPU DSP 948 during the hearing test can be temporarily stored in the RPU DSP 948 for the duration of the test, then deleted upon completion of the test to allow re-use of RPU DSP 948 memory resources during normal operation. The hearing test program may be loaded into the RPU DSP 948</p>

Claims	Prior Art
	<p>through the secondary wireless link 944 or a wired peripheral link 950.</p> <p>Eaton '458 at 6:29–7:6</p> <p>The hearing aid system 102 can communicate with a device 106 facilitated by the short-range network 104. In various embodiments, the device 106 includes a mobile device or a terminal. The term “mobile device” means the inclusion of a digital cellular telephone, a personal digital assistant, a personal communication and information device. Furthermore, the term “mobile device” means the inclusion of a handheld device. The term “terminal” means the inclusion of a data terminal. The term “personal digital assistant” means the inclusion of a portable personal-computing device. The term “personal communication and information device” means the inclusion of a device capable of mobile communication as well as being capable of functioning as a personal digital assistant.</p> <p>Eaton '458 at 11:6–15</p> <p>Figure 6 is a pictorial diagram illustrating a system in accordance with one embodiment. Figure 6 contains elements similar to those discussed in Figure 1. For these elements, the aforementioned discussion related to them is incorporated here in full. The system 600 includes hearing aid systems 602_A and 602_B. The hearing aid system 602_A comprises a programming module 602_{A1} and a hearing aid device 602_{A0}. The hearing aid system 602_B comprises the hearing aid device 602_{B0}. Whereas a mobile device 606 communicates with the hearing aid device 602_{A0} of the hearing aid system 602_A through the programming module 602_{A1}, the mobile device 606 communicates directly with the hearing aid device 602_{B0} of the hearing aid system 602_B.</p> <p>Eaton '458 at 11:16–28</p> <p>For illustrative purposes only, suppose a patient is being fitted with the hearing aid system 602_B. During the fitting process, a piece of software may be executed on the mobile device 606 to interact with the patient wearing the hearing aid system 602_B. Such interaction includes sending audio test signals from the software executing on the mobile device 606 to the hearing aid system 602_B. Such software includes a user interface. The aural response from the patient is either automatically sent back to the testing software by the hearing aid system 602_B or is recorded manually into the testing software through the mobile device 606. Once sufficient aural responses are collected, in one embodiment, such aural responses are sent to a programming fitting server 616 through the cellular/mobile network 610 and the Internet 614; in another embodiment, such aural responses are formatted to form an audiogram before sending the</p>

Claims	Prior Art
	<p>information to the programming fitting server 616.</p> <p>Eaton '458 at 11:28–30 The term audiogram means the inclusion of profiling from the aural responses so as to obtain the extent of the hearing loss of the patient.</p> <p>A POSITA would have been motivated to combine the teachings of the DHA system with Anderson or Eaton, to meet the claim limitations with a reasonable expectation of success for the reasons described above with respect to claims [1b] and 3. In addition, a POSITA would be motivated to combine Anderson and Eaton because the DHA system and Anderson discloses a hearing test for determining signal enhancement parameters for a hearing aid, and Eaton discloses a specific technique for performing such a test using interactive software in order to fit a user with a hearing aid. Eaton, 11:16-33. DHA system and Anderson disclose that the hearing test can be performed using a hearing test program (i.e., software) stored externally. <i>See</i> Anderson, 27:17-20 Thus, the combined DHA system and Anderson lends itself to combination with Eaton's teaching of a specific software program that could be stored externally to perform the hearing test. Moreover, Eaton's hearing aid fitting software would have been readily combinable with the combined DHA system and Anderson's hearing aid system by the POSITA with a reasonable expectation of success because the hearing aid systems of the DHA system, Anderson and Eaton share multiple similarities. For example, each discloses hearing aid systems for use with a wireless communication system and shares analogous and familiar computer components, such as ear-worn hearing aid devices that can be wirelessly connected to a mobile communication device. Accordingly, combining the teachings of Anderson and Eaton would involve combining known components to yield predictable results.</p> <p>To the extent that the DHA system, alone or in view of Anderson and Eaton, does not disclose any elements of this claim limitation, it would have been obvious to a person having ordinary skill in the art to combine, with the teachings of DHA systems, those elements, as known to them or as disclosed by other prior art identified in the cover pleading.</p>
<p>[20b] wherein the mobile communication device is configured to establish the hearing profile associated with the user via the software program executed by the mobile communication device based on the user interacting with</p>	<p>DHA systems, alone or in view of Anderson, expressly or inherently discloses “wherein the mobile communication device is configured to establish the hearing profile associated with the user via the software program executed by the mobile communication device based on the user interacting with the software program via the mobile communication device.”</p> <p><i>See</i> above at claim [20a].</p>

Claims	Prior Art
<p>the software program via the mobile communication device.</p>	<p>A POSITA would have been motivated to combine the teachings of the DHA system with Anderson to meet the claim limitations with a reasonable expectation of success for the reasons described above with respect to claim [1b] and [20a].</p> <p>To the extent that the DHA system, alone or in view of Anderson, does not disclose any elements of this claim limitation, it would have been obvious to a person having ordinary skill in the art to combine, with the teachings of DHA systems, those elements, as known to them or as disclosed by other prior art identified in the cover pleading.</p>
<p>21. A personal audio system as recited in claim 20, wherein the software program at least interactively guides the user to establish the hearing profile for the user.</p>	<p>DHA systems, alone or in view of Anderson and Eaton, expressly or inherently discloses “[a] personal audio system as recited in claim 20, wherein the software program at least interactively guides the user to establish the hearing profile for the user.”</p> <p><i>See above at claim [20a].</i></p> <p>A POSITA would have been motivated to combine the teachings of the DHA system with Anderson to meet the claim limitations with a reasonable expectation of success for the reasons described above with respect to claim [1b] and [20a].</p> <p>To the extent that the DHA system, alone or in view of Anderson and Eaton, does not disclose any elements of this claim limitation, it would have been obvious to a person having ordinary skill in the art to combine, with the teachings of DHA systems, those elements, as known to them or as disclosed by other prior art identified in the cover pleading.</p>
<p>22. A personal audio system as recited in claim 21, wherein the software program at least interactively guides the user to establish the hearing profile for the user by at least causing various audio tones to be output to the user via at least one of the earbuds, prompting the user to provide a hearing input with respect to one or more of the various audio tones being output, and receiving the hearing user input pertaining to</p>	<p>DHA systems, alone or in view of Anderson and Eaton, expressly or inherently discloses “[a] personal audio system as recited in claim 21, wherein the software program at least interactively guides the user to establish the hearing profile for the user by at least causing various audio tones to be output to the user via at least one of the earbuds, prompting the user to provide a hearing input with respect to one or more of the various audio tones being output, and receiving the hearing user input pertaining to one or more of the various audio tones being output.”</p> <p><i>See above at claim 3.</i></p> <p>A POSITA would have been motivated to combine the teachings of the DHA system with Anderson to meet the claim limitations with a reasonable expectation of success for the reasons described above with respect to claim [1b] and [20a].</p> <p>To the extent that the DHA system, alone or in view of Anderson and Eaton, does not disclose any elements of this claim limitation, it would have been obvious to a person having ordinary skill in the art to combine,</p>

Claims	Prior Art
one or more of the various audio tones being output.	with the teachings of DHA systems, those elements, as known to them or as disclosed by other prior art identified in the cover pleading.
<p>[23a] A personal audio system as recited in claim 1, wherein at least one of the earbuds in the pair of earbuds comprises a radiation detector configured to capture radiation data; and</p>	<p>DHA systems, alone or in view of Anderson and Abreu or Irvin '332 or Luo '446, expressly or inherently discloses “[a] personal audio system as recited in claim 1, wherein at least one of the earbuds in the pair of earbuds comprises a radiation detector configured to capture radiation data.”</p> <p>Anderson '783 at 24:15–26</p> <p>The invention provides useful unsolicited information (i.e., information automatically provided to the user by the system, not as a result of a user request for information) to the user of the system in a manner that can be made imperceptible to a casual observer if desired. Information about the status of the hearing aid system itself, such as a verbal warning that an earpiece or RPU battery voltage is low, is provided to the user as detailed earlier. Automatic measurements of heart pulse rate, body temperature and walking distance traveled can be taken by sensors 950 in the RPU 16 or earpiece 10 and stored in the RPU DSP 948.</p> <p>Abreu '960 at 6:28–38</p> <p>Head-mounted structures are structures mounted on the head or neck for positioning sensors on the end of the tunnel and include head bands with accessories that are adjacent to the tunnel, visors, helmets, headphone, structures wrapping around the ear and the like. For the purpose of this description TempAlert is referred herein as a system that measures temperature in the BTT area and has means to report the measured value and that can incorporate alarm means that are activated when certain levels are reached. Support structures yet include any article that has sensing means in which said sensing means are positioned at the end of the tunnel.</p> <p>Abreu '960 at 10:49–55</p> <p>The apparatus of the invention include a temperature sensor located in the medial canthal pad. For temperature measurement the sensing system is located on a skin area that includes the medial canthal corner of the eye and upper eyelid. The sensor in the medial canthal pad is preferably positioned on the skin adjacent to the medial canthal area (medial corner of the eye).</p> <p>Abreu '960 at 32:29–45</p> <p>In accordance with this exemplary embodiment for temperature measurement, the thermal energy emitted by the BTT is sensed by the temperature sensor 102 such as a miniature thermistor which produces a signal representing the thermal energy sensed. The signal is then converted to digital information and</p>

Claims	Prior Art
	<p>processed by processor 106 using standard processing for determining the temperature.</p> <p>Irvin '332 at 6:3–33</p> <p>FIG. 4C illustrates the power control circuit 210 coupled with a contact switch functioning as the sensor 220. The contact switch includes a first and second contact connected to the power control circuit 210, and connecting member 222 actuated when the user places the wireless headset 100 on their body. When actuated, the connecting member 222 places the first and second contacts of sensor 220 in electrical connection. Ideally, the sensor 220, when configured as a contact sensor, is placed in the earpiece 130, in the retaining member 140, or in the retaining pad 150 shown in FIG. 1B. Preferably, when sensor 220 is configured as a contact-based proximity sensor, the contacting element of sensor 220 is sealed in a manner that minimizes contaminant ingress. In some embodiments, the contact sensor may be actuated by the user resiliently expanding the retaining member 140 or other element of the wireless headset 100 for mounting upon their body.</p> <p>FIG. 4D illustrates the power control circuit 210 coupled with a proximity sensor 220 not based on mechanical actuation. The proximity sensor 220 may be an active or passive IR sensor, may be a resistive or other type of thermal sensor, or may be an acoustic sensor. As a proximity sensor, sensor 220 provides indication that a user has placed the wireless headset 100 upon their body. As such, any of the foregoing or other presence sensing techniques may be used to implement sensor 220. Preferably, when configured as a proximity sensor, sensor 220 requires minimal power and provides a sensor output signal based on proximity to the user's body. Depending upon the sensing technology chosen, the control circuit 210 may be required to periodically stimulate sensor 220 to generate a detection response.</p> <p>Luo '446 at [0036]</p> <p>The invention generally relies on shining light having a particular wavelength on human skin as well as the reflectance properties of human skin. The surface reflection of the human skin occurs at the surface of the epidermis and is approximately independent of the lighting wavelength and independent of human race. In the infrared (IR) range the wavelengths that can be used are greater than 800 nano-meters and the surface reflectance of human skin is more or less constant and close to 50%. The same hold true for orange and red light in the visible light range, namely wavelengths between approximately 600 and 800 nano-meters can be used, although the surface reflectance for darker skin is reduced for shorter wavelengths. Surface reflectance for dark skin is approximately 25% at 700 nano-meters. In the embodiments shown herein, the energy of choice in the sensor unit 18 is preferably infrared (IR) energy. However, long wavelength visible light energy can also be used as discussed below.</p>

Claims	Prior Art
	<p>A POSITA would have been motivated to combine the teachings of the DHA system with Anderson and Abreu or Luo or Irvin to meet the claim limitations with a reasonable expectation of success for the reasons described above with respect to claims [1a]. In addition, a POSITA would have been motivated to combine the teachings of the DHA Systems, such as the Phonak system, with Anderson and Abreu or Luo or Irvin to meet the claim limitations with a reasonable expectation of success. For example, both the DHA system and Abreu disclose head-worn devices. The Phonak system, for example, discloses that its device can include a sensor (a microphone), so it would have been obvious to the POSITA to incorporate other types of sensors into the device as well. Given that Abreu discloses a head-worn device with various sensors for measuring different types of data, the POSITA would have looked to Abreu or Luo or Irvin for different types of sensors that could be incorporated into the Phonak system. Moreover, a POSITA would have understood that sending a signal if the earpiece is no longer being worn by a user would necessarily include processing the data from the sensor to make that determination. Further, a POSITA would have understood that the sensor would fail to capture signals in a proper manner if the earbud were not being worn. Accordingly, the modified DHA system could determine whether the Phonak system's earbuds are being worn based on at least the sensor data. Further, a POSITA would have been motivated to modify the DHA system to include this feature to conserve battery use, as if the earbud is not being worn, the device can switch to a low power mode and minimize power usage. Moreover, the devices of the DHA system, Anderson '783 and Abreu or Luo or Irvin share various structural and functional similarities that render the inventions readily combinable. Accordingly, combining the teachings of the DHA system with Anderson '783 and Abreu or Luo or Irvin would involve combining known components to yield predictable results.</p> <p>To the extent that the DHA system, alone or in view of Anderson and Abreu or Luo or Irvin, does not disclose any elements of this claim limitation, it would have been obvious to a person having ordinary skill in the art to combine, with the teachings of DHA systems, those elements, as known to them or as disclosed by other prior art identified in the cover pleading.</p>
<p>[23b] first electronic circuitry electrically connected to the radiation detector, the first electronic circuitry configured to process at least the radiation data captured by the radiation detector,</p>	<p>DHA systems, alone or in view of Anderson and Abreu or Irvin '332 or Luo '446, expressly or inherently discloses "first electronic circuitry electrically connected to the radiation detector, the first electronic circuitry configured to process at least the radiation data captured by the radiation detector."</p> <p>Anderson '783 at 12:47–13:22</p> <p>Remote processor unit interrogator and reply processor: The RPU, shown in FIG. 9, may use any of a wide variety of components well known in the art. The DSP 948 may contain several integrated circuits, e.g., the Motorola DSP56L002 single-chip digital signal processor, a Motorola DSP56200 cascadable-adaptive digital filter chip and supporting memory devices well known in the art. Miniature components</p>

Claims	Prior Art
	<p>such as dielectric resonance oscillators are available for the microwave oscillator 908, and the DPSK replies can be decoded using surface acoustic wave (SAW) or glass bulk acoustic wave (BAW) delay lines 936 in a standard delay-and-multiply demodulator circuit 934. An alternative demodulator that achieves superior noise immunity can be implemented using a phase-locked loop in a squaring loop configuration for carrier recovery instead of delay-and-multiply circuitry 934. Decoding of a DPSK waveform using such coherent detection techniques is known in the art as differentially encoded coherent phase shift keying, or DECPSK. Amplifiers using low noise pseudomorphic high electron mobility transistors (e.g., the HP ATF-36077) may be used as needed (e.g., between the circulator 906 and routing switch 922), and such amplifiers are not individually identified in FIG. 9. All switches shown in FIG. 9 are RF switches that can be implemented using, e.g., PIN diodes. Switch control connections from the DSP 948 to the RF switches are shown in FIG. 9 for clarity, and all switches remain set as shown until changed as noted. The system shown in FIG. 9 demodulates a received reply when the RF reply travels through a selected antenna 900, antenna selection switch 902, circulator 906, routing switch 922 (set to the position opposite that shown in FIG. 9 during the reply interval by DSP 948 via control signal 976), bandpass filter 928 and IF mixer 930. The IF mixer 930 downconverts the reply by mixing it with a small amount of RF power from the 5760 MHz oscillator 908 via switch 910 (which is set to the position opposite that shown in FIG. 9 during the reply interval by DSP 948 via control signal 972) and power splitter 920. The resulting DPSK reply has a center frequency of 60 MHz, a bandwidth of 100 MHz, and is passed through the bandpass filter 932 to the DPSK demodulator 934 which, in this example, consists of a 50 nsec delay line 936 and multiplier 938. The resulting bits are stored in the reply data register 940 for processing by the DSP 948.</p> <p>Anderson '783 at 24:15–26</p> <p>The invention provides useful unsolicited information (i.e., information automatically provided to the user by the system, not as a result of a user request for information) to the user of the system in a manner that can be made imperceptible to a casual observer if desired. Information about the status of the hearing aid system itself, such as a verbal warning that an earpiece or RPU battery voltage is low, is provided to the user as detailed earlier. Automatic measurements of heart pulse rate, body temperature and walking distance traveled can be taken by sensors 950 in the RPU 16 or earpiece 10 and stored in the RPU DSP 948.</p> <p>Anderson '783 at 24:26–32</p> <p>Sensor measurements can be reported to the user by a voice synthesized in the RPU DSP 948 using any of a variety of synthesis techniques well known in the art, including playback of pre-recorded digitized</p>

Claims	Prior Art
	<p>speech. The voice signal produced within the RPU DSP 948 is then added to the processed audio signals already being sent from the RPU to the earpiece.</p> <p>Abreu '960 at 6:28–38 Head-mounted structures are structures mounted on the head or neck for positioning sensors on the end of the tunnel and include head bands with accessories that are adjacent to the tunnel, visors, helmets, headphone, structures wrapping around the ear and the like. For the purpose of this description TempAlert is referred herein as a system that measures temperature in the BTT area and has means to report the measured value and that can incorporate alarm means that are activated when certain levels are reached. Support structures yet include any article that has sensing means in which said sensing means are positioned at the end of the tunnel.</p> <p>Abreu '960 at 10:49–55 The apparatus of the invention include a temperature sensor located in the medial canthal pad. For temperature measurement the sensing system is located on a skin area that includes the medial canthal corner of the eye and upper eyelid. The sensor in the medial canthal pad is preferably positioned on the skin adjacent to the medial canthal area (medial corner of the eye).</p> <p>Abreu '960 at 15:57–60 The eyeglasses can be incorporated with a detecting system to send a signal if the eyeglasses were removed or if the temperature sensor is not capturing signals in a proper manner.</p> <p>Abreu '960 at 18:56–59 The present invention preferably includes all functions in a miniature semiconductor chip, which as an integrated circuit, incorporates sensor, processing and transmitting units and control circuits.</p> <p>Abreu '960 at 32:29–45 In accordance with this exemplary embodiment for temperature measurement, the thermal energy emitted by the BTT is sensed by the temperature sensor 102 such as a miniature thermistor which produces a signal representing the thermal energy sensed. The signal is then converted to digital information and processed by processor 106 using standard processing for determining the temperature.</p>

Claims	Prior Art
	<p>Irvin '332 at 6:3–33</p> <p>FIG. 4C illustrates the power control circuit 210 coupled with a contact switch functioning as the sensor 220. The contact switch includes a first and second contact connected to the power control circuit 210, and connecting member 222 actuated when the user places the wireless headset 100 on their body. When actuated, the connecting member 222 places the first and second contacts of sensor 220 in electrical connection. Ideally, the sensor 220, when configured as a contact sensor, is placed in the earpiece 130, in the retaining member 140, or in the retaining pad 150 shown in FIG. 1B. Preferably, when sensor 220 is configured as a contact-based proximity sensor, the contacting element of sensor 220 is sealed in a manner that minimizes contaminant ingress. In some embodiments, the contact sensor may be actuated by the user resiliently expanding the retaining member 140 or other element of the wireless headset 100 for mounting upon their body.</p> <p>FIG. 4D illustrates the power control circuit 210 coupled with a proximity sensor 220 not based on mechanical actuation. The proximity sensor 220 may be an active or passive IR sensor, may be a resistive or other type of thermal sensor, or may be an acoustic sensor. As a proximity sensor, sensor 220 provides indication that a user has placed the wireless headset 100 upon their body. As such, any of the foregoing or other presence sensing techniques may be used to implement sensor 220. Preferably, when configured as a proximity sensor, sensor 220 requires minimal power and provides a sensor output signal based on proximity to the user's body. Depending upon the sensing technology chosen, the control circuit 210 may be required to periodically stimulate sensor 220 to generate a detection response.</p> <p>Luo '446 at [0034]</p> <p>In general, the location sensor module 18 and the hearing aid module 24 communicate via a bi-directional information signal 42. The hearing aid module 24 polls the location sensor module 18 via the bi-directional signal 42, on preferably a periodic basis, to determine whether the hearing aid 10 is in the ear of the hearing aid user. In this case, the bi-directional signal 42 acts as a polling signal. In response to polling done via the bi-directional information signal 42, the location sensor module 18 probes the outer environment of the hearing aid 10 and returns location information via the bi-directional location signal 42. In this case, the bi-directional signal 42 acts as a location information signal. If the location information indicates that the hearing aid 10 is in the ear of the hearing aid user (i.e. the in-the-ear case), the hearing aid 10 operates in full-function mode. If the location information indicates that the hearing aid 10 is not in the ear of the hearing aid user (i.e. the out-of-the-ear case), the hearing aid 10 operates in sleep mode. In sleep mode, the hearing aid 10 can essentially be considered to be off.</p> <p>Luo '446 at [0035]</p>

Claims	Prior Art
	<p>There are several general scenarios for the location of the hearing aid. In the first scenario, the hearing aid 10 is not in the ear of the user and the hearing aid 10 is in sleep mode. In this case, the location sensor module 18 is polled and the location information indicates that the hearing aid 10 is not in the ear of the hearing aid user and the hearing aid 10 continues to operate in sleep mode. In the second scenario, the hearing aid 10 has just been inserted into the ear of the hearing aid user and the hearing aid 10 was previously in sleep mode. In this case, the location sensor module 18 is polled and the location information indicates that the hearing aid 10 is in the ear of the hearing aid user. The hearing aid 10 then moves into full-function mode. In the third scenario, the hearing aid 10 is in full-function mode and is being taken out of the ear of the hearing aid user. The location sensor module 18 is polled and the location information indicates that the hearing aid 10 is no longer in the ear of the hearing aid user. The hearing aid 10 then moves into sleep mode.</p> <p>A POSITA would have been motivated to combine the teachings of the DHA system with Anderson and Abreu or Luo or Irvin to meet the claim limitations with a reasonable expectation of success for the reasons described above with respect to claim [1a] and [23a].</p> <p>To the extent that the DHA system, alone or in view of Anderson and Abreu or Luo or Irvin, does not disclose any elements of this claim limitation, it would have been obvious to a person having ordinary skill in the art to combine, with the teachings of DHA systems, those elements, as known to them or as disclosed by other prior art identified in the cover pleading.</p>
<p>[23c] wherein the first electronic circuitry is configured to determine whether the earbud is being worn by a user based on at least the processed radiation data.</p>	<p>DHA systems, alone or in view of Anderson and Abreu or Luo or Irvin, expressly or inherently discloses “wherein the first electronic circuitry is configured to determine whether the earbud is being worn by a user based on at least the processed radiation data.”</p> <p>Abreu ’960 at 15:57–60 The eyeglasses can be incorporated with a detecting system to send a signal if the eyeglasses were removed or if the temperature sensor is not capturing signals in a proper manner.</p> <p>Abreu ’960 at 18:56–59 The present invention preferably includes all functions in a miniature semiconductor chip, which as an integrated circuit, incorporates sensor, processing and transmitting units and control circuits.</p> <p>Abreu ’960 at 44:61–67 The present invention includes means for directing preferably near-infrared energy into the surface of the skin at the end of the BTE tunnel, means for analyzing and converting the reflectance or back scattered</p>

Claims	Prior Art
	<p>spectrum into the concentration of the substance measured and support means for positioning the light source and detector means adjacent to the surface of the skin at the BTE tunnel site.</p> <p>Irvin '332 at 2:20–34</p> <p>The wireless headset serves as a communication accessory to a base unit, such as a cellular telephone. When not in use, the wireless headset reverts to an inactive state and draws essentially no power. A power control circuit imparts automatic enabling capability to the wireless headset. Preferably, the power control circuit employs a sensor responsive to a natural stimulus associated with a person preparing to use the wireless headset. The wireless headset may sense motion, physical orientation, or user proximity. Sensing techniques include but are not limited to micro-electro-mechanical-systems (MEMS), mercury switches, attitude switches, and IR or acoustic proximity sensors. In some embodiments, the wireless headset transitions from the inactive state to an active state in response to a user picking it up or mounting the wireless headset upon their head.</p> <p>Irvin '332 at 6:3–33</p> <p>FIG. 4C illustrates the power control circuit 210 coupled with a contact switch functioning as the sensor 220. The contact switch includes a first and second contact connected to the power control circuit 210, and connecting member 222 actuated when the user places the wireless headset 100 on their body. When actuated, the connecting member 222 places the first and second contacts of sensor 220 in electrical connection. Ideally, the sensor 220, when configured as a contact sensor, is placed in the earpiece 130, in the retaining member 140, or in the retaining pad 150 shown in FIG. 1B. Preferably, when sensor 220 is configured as a contact-based proximity sensor, the contacting element of sensor 220 is sealed in a manner that minimizes contaminant ingress. In some embodiments, the contact sensor may be actuated by the user resiliently expanding the retaining member 140 or other element of the wireless headset 100 for mounting upon their body.</p> <p>FIG. 4D illustrates the power control circuit 210 coupled with a proximity sensor 220 not based on mechanical actuation. The proximity sensor 220 may be an active or passive IR sensor, may be a resistive or other type of thermal sensor, or may be an acoustic sensor. As a proximity sensor, sensor 220 provides indication that a user has placed the wireless headset 100 upon their body. As such, any of the foregoing or other presence sensing techniques may be used to implement sensor 220. Preferably, when configured as a proximity sensor, sensor 220 requires minimal power and provides a sensor output signal based on proximity to the user's body. Depending upon the sensing technology chosen, the control circuit 210 may be required to periodically stimulate sensor 220 to generate a detection response.</p>

Claims	Prior Art
	<p data-bbox="583 175 861 207">Irvin '332 at 8:11–22</p> <p data-bbox="632 215 1992 467">FIGS. 6 and 7 illustrate operation of the wireless headset 100. FIG. 6 is a state diagram for one exemplary embodiment of the present invention. When the wireless headset 100 is not in use, it automatically enters an inactive state (610). As explained, this inactive state may constitute disconnecting the communications circuit 230 from the power source 240, or may simply represent placing the communications circuit 230 in a non-functional, low-power condition. Absent user stimulus, such as motion of the wireless headset 100 or placement of the wireless headset 100 upon the user's body, the wireless headset 100 remains in the inactive state (610).</p> <p data-bbox="583 524 827 557">Luo '446 at [0034]</p> <p data-bbox="632 565 1992 963">In general, the location sensor module 18 and the hearing aid module 24 communicate via a bi-directional information signal 42. The hearing aid module 24 polls the location sensor module 18 via the bi-directional signal 42, on preferably a periodic basis, to determine whether the hearing aid 10 is in the ear of the hearing aid user. In this case, the bi-directional signal 42 acts as a polling signal. In response to polling done via the bi-directional information signal 42, the location sensor module 18 probes the outer environment of the hearing aid 10 and returns location information via the bi-directional location signal 42. In this case, the bi-directional signal 42 acts as a location information signal. If the location information indicates that the hearing aid 10 is in the ear of the hearing aid user (i.e. the in-the-ear case), the hearing aid 10 operates in full-function mode. If the location information indicates that the hearing aid 10 is not in the ear of the hearing aid user (i.e. the out-of-the-ear case), the hearing aid 10 operates in sleep mode. In sleep mode, the hearing aid 10 can essentially be considered to be off.</p> <p data-bbox="583 1019 827 1052">Luo '446 at [0035]</p> <p data-bbox="632 1060 1992 1409">There are several general scenarios for the location of the hearing aid. In the first scenario, the hearing aid 10 is not in the ear of the user and the hearing aid 10 and is in sleep mode. In this case, the location sensor module 18 is polled and the location information indicates that the hearing aid 10 is not in the ear of the hearing aid user and the hearing aid 10 continues to operate in sleep mode. In the second scenario, the hearing aid 10 has just been inserted into the ear of the hearing aid user and the hearing aid 10 was previously in sleep mode. In this case, the location sensor module 18 is polled and the location information indicates that the hearing aid 10 is in the ear of the hearing aid user. The hearing aid 10 then moves into full-function mode. In the third scenario, the hearing aid 10 is in full-function mode and is being taken out of the ear of the hearing aid user. The location sensor module 18 is polled and the location information indicates that the hearing aid 10 is no longer in the ear of the hearing aid user. The hearing aid</p>

Claims	Prior Art
	<p>10 then moves into sleep mode.</p> <p>Luo '446 at [0036] The invention generally relies on shining light having a particular wavelength on human skin as well as the reflectance properties of human skin. The surface reflection of the human skin occurs at the surface of the epidermis and is approximately independent of the lighting wavelength and independent of human race. In the infrared (IR) range the wavelengths that can be used are greater than 800 nano-meters and the surface reflectance of human skin is more or less constant and close to 50%. The same hold true for orange and red light in the visible light range, namely wavelengths between approximately 600 and 800 nano-meters can be used, although the surface reflectance for darker skin is reduced for shorter wavelengths. Surface reflectance for dark skin is approximately 25% at 700 nano-meters. In the embodiments shown herein, the energy of choice in the sensor unit 18 is preferably infrared (IR) energy. However, long wavelength visible light energy can also be used as discussed below.</p> <p>A POSITA would have been motivated to combine the teachings of the DHA system with Anderson and Abreu or Lu or Irvin to meet the claim limitations with a reasonable expectation of success for the reasons described above with respect to claim [1a] and [23a].</p> <p>To the extent that the DHA system, alone or in view of Anderson and Abreu, does not disclose any elements of this claim limitation, it would have been obvious to a person having ordinary skill in the art to combine, with the teachings of DHA systems, those elements, as known to them or as disclosed by other prior art identified in the cover pleading.</p>
24[pre]. A personal audio system operable to at least generate customized audio output, comprising:	DHA systems, alone or in view of Anderson and/or Enzmann, expressly or inherently discloses “[a] personal audio system operable to at least generate customized audio output.” <i>See</i> above at claim 1[pre]. A POSITA would have been motivated to combine the teachings of the DHA system with Anderson and/or Enzmann to meet the claim limitations with a reasonable expectation of success for the reasons described above with respect to claim 1. To the extent that the DHA system, alone or in view of Anderson, does not disclose any elements of this claim limitation, it would have been obvious to a person having ordinary skill in the art to combine, with the teachings of DHA systems, those elements, as known to them or as disclosed by other prior art identified in the cover pleading.

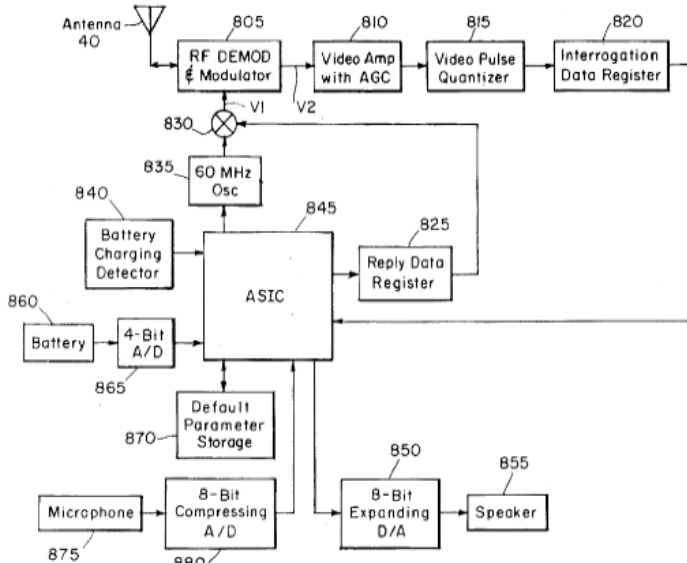
Claims	Prior Art
<p>[24a] a head-worn wearable electronic device including at least a pair of speakers, each of the speakers being configured to be wearable proximate to an ear of a user; and</p>	<p>DHA systems, alone or in view of Anderson, expressly or inherently discloses “a head-worn wearable electronic device including at least a pair of speakers, each of the speakers being configured to be wearable proximate to an ear of a user.”</p> <p><i>See</i> above at claim [1a].</p> <p>Moreover, Anderson ’783 discloses that the earpiece of a hearing aid includes a speaker, and can be worn “completely-in-the-canal” of the ear, as shown in FIG. 1 or behind the ear, or formed from a headset 21. <i>See id.</i> at 3:52-60; 4:15-39; 5:52-55; FIG. 1; FIG. 2. The speaker is thus wearable proximate to the ear of the user.</p> <p>In some embodiments, such as shown in FIG. 2, the earpiece belongs to a head-worn headset 21. <i>See</i> 4:15-39; FIG. 2.</p> <div data-bbox="592 574 1533 1149" data-label="Diagram"> <pre> graph LR subgraph Headset_21 [Headset 21] WM[Wireless Microphone 20] end subgraph Remote_Processor_Unit_23 [Remote Processor Unit 23] FM_XMTR[FM XMTR 27] FM_RCVR[FM RCVR 24] SE[Signal Enhancer 26] DPDT_Switch_25[DPDT Switch 25] SPKR[SPKR] MIC[MIC] FM_XMTR --> SE SE --> SPKR FM_RCVR --> DPDT_Switch_25 DPDT_Switch_25 --> MIC DPDT_Switch_25 --> SPKR end subgraph Handset_28 [Cordless or Cellular Telephone Handset 28] SWL[Secondary Wireless Link 29] end Headset_21 -- F1 --> Remote_Processor_Unit_23 Remote_Processor_Unit_23 -- F2 --> Headset_21 Remote_Processor_Unit_23 <--> Secondary Wireless Link 29 Handset_28 </pre> </div> <p>FIG. 2</p> <p>The headset earpiece 22 shown in FIG. 2 is “functionally similar” to the in-ear earpiece 10 shown in FIG. 1. <i>See</i> 4:33-35. In a disclosed embodiment, “[a]n FM wireless microphone 20 transmitter, e.g. Radio Shack (R) model 33-1076, is mounted on an FM headset 21 receiver, e.g. Radio Shack (R) model 12-103, to form an earpiece 22 which is not necessarily hidden from view.” <i>See</i> 4:33-39. A POSITA would understand such Radio Shack models to include a headset that is “head-worn.” The illustrated headset 21, as well as the</p>

Claims	Prior Art
	<p>referenced commercially available models, includes a headband that rests on the top of the user's head for support.</p> <p>A POSITA would have understood from these disclosures that Anderson's headset 21 can include a pair of speakers, one for each ear, as shown above; with each speaker wearable near (and thus proximate to) its respective ear, <i>i.e.</i>, when the headset is worn.</p> <p>Anderson illustrates headset 21 and earpiece 22 within the context of a "limited-feature embodiment," shown in FIG. 2. <i>See</i> 4:15-29. However, a POSITA would have understood headset 21 and earpiece 22 to be combinable with the "full-featured" embodiment of FIG. 1 and/or DHA Systems. <i>See</i> 3:52-58. This is because Anderson characterizes FIG. 1 and FIG. 2 as showing two preferred embodiments of the same invention; explains that the headset earpiece shown in FIG. 2 is "functionally similar" to the in-ear earpiece 10 of FIG. 1; and further explains that "many other earpiece housing types" can be used in place of the CIC earpiece housing shown in FIG. 1. <i>See</i> 3:52-58, 4:15-35. A POSITA would have understood from these disclosures that the FIG. 1 system and/or DHA Systems could be flexibly implemented using various housings for the earpiece, including headset 21.</p> <p>A POSITA would have been motivated to use the headset instead of the CIC housing or the housing of the DHA Systems to realize advantages in cost, simplicity, and commercial availability described by Anderson. Even using headset 21 in a full-featured application and/or DHA Systems, it would provide certain advantages over the in-ear embodiment. Unlike some CIC housings, a headset does not need to be customized to fit in a user's ear canal. Additionally, there would have been comfort benefits: many users do not like supporting objects with their ears, a problem addressed by replacing a CIC housing or DHA System housing with a headset. And using a larger headset rather than a small in-ear device would make it easier and less expensive to include the processing and communication capabilities of the full-featured application, and would reduce the likelihood of misplacing the device. Anderson provides specific guidance for adapting the headset using "techniques [that] are well known in the art," and notes that schematics and service manuals for the headset were readily available, demonstrating that the POSITA would have had a reasonable expectation of success in making the combination. <i>See</i> 4:26-5:60.</p> <p>Because the POSITA would have combined headset 21 and its earpiece 22 with the full-featured FIG. 1 embodiment, throughout these invalidity contentions for the '693 and '355 patents, unless otherwise noted, references to the earpiece 10 shown for FIG. 1 should be understood to apply equally to earpiece 22 and headset 21.</p>

Claims	Prior Art
	<p>The “head-worn wearable electronic device” is “eyewear or headwear” and thus encompasses an article worn over the user’s head in the manner shown in Fig. 2. In addition, it would have been obvious to a POSITA to make the in-ear earpiece of Anderson and/or the DHA Systems into a head–worn device, e.g., to disguise the earpiece and/or to reduce the number of components a user needs to wear for different purposes, and in particular in view of the well known nature of head–worn wearable electronic devices at the time of the invention.</p> <p>In addition, as noted above, Anderson ’783 discloses that “it will be clear to those skilled in the art that a binaural (stereo) system can be constructed based on the same principles as the monaural system.” <i>See id.</i> at 4:30-33. A POSITA would thus have understood the disclosure to encompass a pair of earbuds with a pair of speakers, one for each ear.</p> <p>A POSITA would have been motivated to combine the teachings of the DHA system with Anderson to meet the claim limitations with a reasonable expectation of success for the reasons described above with respect to claim 1.</p> <p>To the extent that the DHA system, alone or in view of Anderson, does not disclose any elements of this claim limitation, it would have been obvious to a person having ordinary skill in the art to combine, with the teachings of DHA systems, those elements, as known to them or as disclosed by other prior art identified in the cover pleading.</p>
<p>[24b] a mobile communication device at least for supporting audio output, the mobile communication device configured to store at least a hearing profile associated with the user, the stored hearing profile associated with the user at least characterizes hearing impairment regarding the user's ability to hear,</p>	<p>DHA systems, alone or in view of Anderson and/or Enzmann, expressly or inherently discloses “a mobile communication device at least for supporting audio output, the mobile communication device configured to store at least a hearing profile associated with the user, the stored hearing profile associated with the user at least characterizes hearing impairment regarding the user's ability to hear.”</p> <p><i>See</i> above at claim [1b].</p> <p>A POSITA would have been motivated to combine the teachings of the DHA system with Anderson and/or Enzmann to meet the claim limitations with a reasonable expectation of success for the reasons described above with respect to claim 1.</p> <p>To the extent that the DHA system, alone or in view of Anderson, does not disclose any elements of this claim limitation, it would have been obvious to a person having ordinary skill in the art to combine, with the teachings of DHA systems, those elements, as known to them or as disclosed by other prior art identified in the cover pleading.</p>
<p>[24c] the mobile communication device configured to produce</p>	<p>DHA systems, alone or in view of Anderson, expressly or inherently discloses “the mobile communication device configured to produce audio output signals customized for the user based on at least the stored hearing</p>

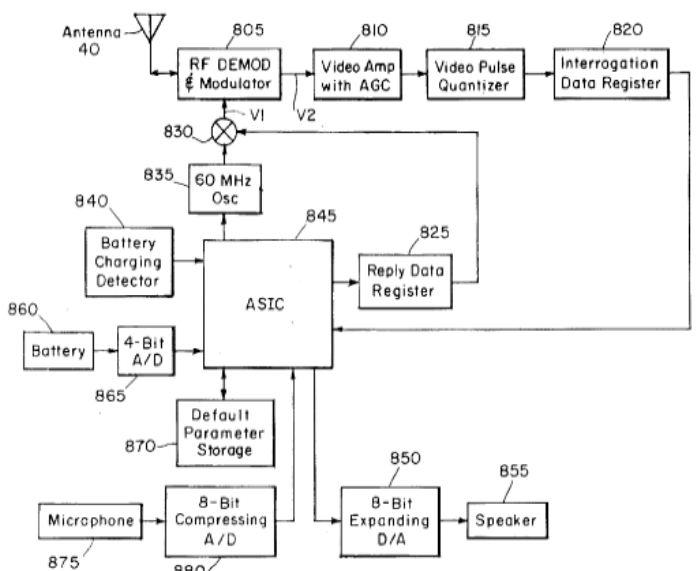
Claims	Prior Art
<p>audio output signals customized for the user based on at least the stored hearing profile associated with the user such that the customized audio output signals at least compensate for the hearing impairment regarding the user's ability to hear,</p>	<p>profile associated with the user such that the customized audio output signals at least compensate for the hearing impairment regarding the user's ability to hear.”</p> <p><i>See above at claim [1c].</i></p> <p>A POSITA would have been motivated to combine the teachings of the DHA system with Anderson to meet the claim limitations with a reasonable expectation of success for the reasons described above with respect to claim 1.</p> <p>To the extent that the DHA system, alone or in view of Anderson, does not disclose any elements of this claim limitation, it would have been obvious to a person having ordinary skill in the art to combine, with the teachings of DHA systems, those elements, as known to them or as disclosed by other prior art identified in the cover pleading.</p>
<p>[24d] the mobile communication device configured to produce the customized audio output signals using at least signal amplification,</p>	<p>DHA systems, alone or in view of Anderson, expressly or inherently discloses “the mobile communication device configured to produce the customized audio output signals using at least signal amplification.”</p> <p><i>See above at claim [1d].</i></p> <p>A POSITA would have been motivated to combine the teachings of the DHA system with Anderson to meet the claim limitations with a reasonable expectation of success for the reasons described above with respect to claim 1.</p> <p>To the extent that the DHA system, alone or in view of Anderson, does not disclose any elements of this claim limitation, it would have been obvious to a person having ordinary skill in the art to combine, with the teachings of DHA systems, those elements, as known to them or as disclosed by other prior art identified in the cover pleading.</p>
<p>[24e] and the mobile communication device configured to transmit the customized audio output signals to the pair of speakers,</p>	<p>DHA systems, alone or in view of Anderson, expressly or inherently discloses “and the mobile communication device configured to transmit the customized audio output signals to the pair of speakers.”</p> <p><i>See above at claim [1e].</i></p> <p>A POSITA would have been motivated to combine the teachings of the DHA system with Anderson to meet the claim limitations with a reasonable expectation of success for the reasons described above with respect to claim 1.</p> <p>To the extent that the DHA system, alone or in view of Anderson, does not disclose any elements of this claim limitation, it would have been obvious to a person having ordinary skill in the art to combine, with the teachings of DHA systems, those elements, as known to them or as disclosed by other prior art identified in the cover pleading.</p>

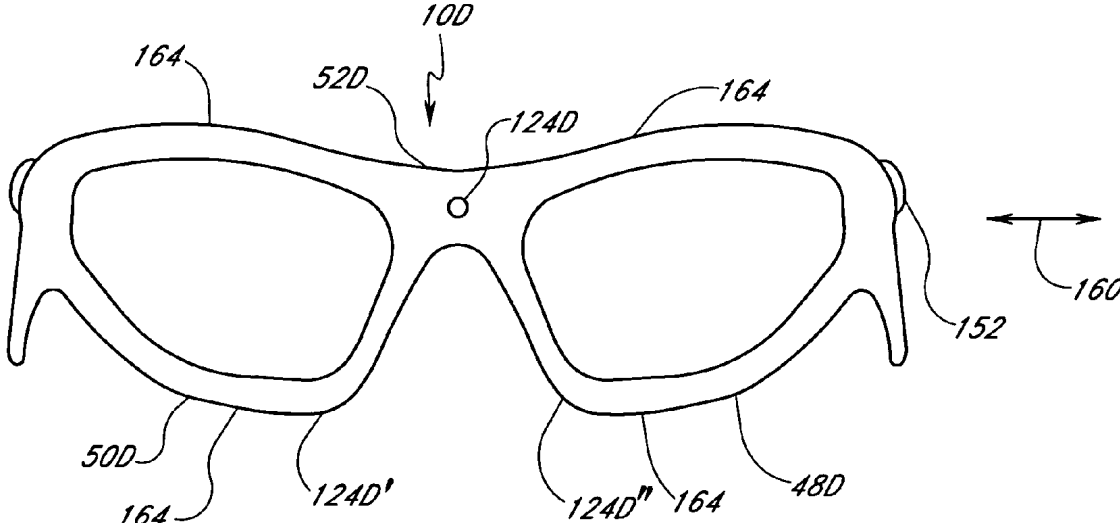
Claims	Prior Art
<p>[24f] wherein each of the speakers is configured to receive at least a portion of the customized audio output signals transmitted by the mobile communication device, and</p>	<p>DHA systems, alone or in view of Anderson, expressly or inherently discloses “wherein each of the speakers is configured to receive at least a portion of the customized audio output signals transmitted by the mobile communication device.”</p> <p><i>See</i> above at claim [1f].</p> <p>A POSITA would have been motivated to combine the teachings of the DHA system with Anderson to meet the claim limitations with a reasonable expectation of success for the reasons described above with respect to claim 1.</p> <p>To the extent that the DHA system, alone or in view of Anderson, does not disclose any elements of this claim limitation, it would have been obvious to a person having ordinary skill in the art to combine, with the teachings of DHA systems, those elements, as known to them or as disclosed by other prior art identified in the cover pleading.</p>
<p>[24g] wherein each of the speakers is configured to generate at least a portion of the customized audio output for the user in accordance with at least the portion of the customized audio output signals being received.</p>	<p>DHA systems, alone or in view of Anderson, expressly or inherently discloses “wherein each of the speakers is configured to generate at least a portion of the customized audio output for the user in accordance with at least the portion of the customized audio output signals being received.”</p> <p><i>See</i> above at claim [1g].</p> <p>A POSITA would have been motivated to combine the teachings of the DHA system with Anderson to meet the claim limitations with a reasonable expectation of success for the reasons described above with respect to claim 1.</p> <p>To the extent that the DHA system, alone or in view of Anderson, does not disclose any elements of this claim limitation, it would have been obvious to a person having ordinary skill in the art to combine, with the teachings of DHA systems, those elements, as known to them or as disclosed by other prior art identified in the cover pleading.</p>
<p>[25a] A personal audio system as recited in claim 24, wherein the head-worn wearable electronic device includes at least a rechargeable battery, and</p>	<p>DHA systems, alone or in view of Anderson, expressly or inherently discloses “[a] personal audio system as recited in claim 24, wherein the head-worn wearable electronic device includes at least a rechargeable battery.”</p> <p>Phonak MicroLink MLxS Product Information at Page. 1</p> <p><small>* Powered by the hearing instrument battery</small></p>

Claims	Prior Art
	<p data-bbox="583 175 898 207">Anderson '783 at FIG. 8</p>  <p data-bbox="945 795 1029 828">FIG. 8</p> <p data-bbox="583 844 919 876">Anderson '783 at 3:40–41</p> <p data-bbox="630 885 1963 950">FIG. 8 is a block diagram showing details of an earpiece for a full-featured preferred embodiment of the invention.</p> <p data-bbox="583 1006 919 1039">Anderson '783 at 6:55–57</p> <p data-bbox="630 1047 1927 1112">Earpiece power consumption is preferably minimized to extend battery life, and use of a rechargeable battery is desirable.</p> <p data-bbox="583 1169 913 1201">Jannard '526 at 20:57–64</p> <p data-bbox="630 1209 1984 1274">The interface device 110 includes a power source 112, a transceiver 114, an interface 116, and an antenna 118.</p> <p data-bbox="630 1282 1953 1347">The power source 112 can be in the form of disposable or rechargeable batteries. Optionally, the power source 112 can be in the form of solar panels and a power regulator.</p> <p data-bbox="630 1356 1848 1421">The transceiver 114 can be in the form of a digital wireless transceiver for one-way or two-way communication.</p>

Claims	Prior Art
	<p>Jannard '526 at 24:37–55</p> <p>Optionally, the power source 112 can comprise portions of the ear stems 54D, 56D which have been formed into batteries. For example, the rear portions 160, 162 of the ear stems 54D, 56D, respectively, can be in the form of custom made batteries, either disposable or rechargeable. Preferably, the rear portions 160, 162 are removable from the forward portions of the ear stems 54D, 56D. This provides a particular advantage in terms of balance. As noted above, imbalanced loads on the head can cause muscular pain and/or headaches. In particular, excessive pressure on the nose can cause severe headaches. Additionally, batteries can have a significantly higher mass density than plastic and lightweight metals, such as sintered titanium or magnesium. Thus, by constructing the rearward portions 160, 162 of the ear stems 54D, 56D of batteries, the weight of these batteries can improve forward-rearward balance of the audio device 10D in that the weight of the interface device 110 can be offset by the batteries. In another embodiment, the ear stems 54D, 56D can define a housing for removable batteries.</p> <p>Jannard '526 at 24:56–25:3</p> <p>The audio device 10D can also include power contacts 164 for recharging any rechargeable batteries connected thereto. For example, the power contacts 164 can be disposed on a lower edge of the orbitals 48D, 50D. Thus, with an appropriate recharging cradle (not shown), the audio device 10D can be laid on the cradle, thereby making contact between the power contacts 164 and corresponding contacts in the cradle (not shown). Alternatively, power contacts can be provided in numerous other locations as desired. For example, the power contacts 164 can be disposed at the ends of the ear stems 54D, 56D. A corresponding cradle can include two vertically oriented holes into which the ear stems are inserted for recharging. In this configuration, the lens within the orbitals 48D, 50D would face directly upwardly.</p> <p>To the extent that The Phonak Digital Hearing Aid system and/or other DHA Systems do not disclose any elements of this claim limitation, it would have been obvious to a person having ordinary skill in the art to combine, with the teachings of The Phonak Digital Hearing Aid system and/or other DHA Systems, those elements, as known to them or as disclosed by other prior art identified in the cover pleading, including Anderson or Jannard. As explained in in reference to limitation [1b], a POSITA would have been motivated to combine the DHA system with the features of secondary references, including Anderson, to meet the claim limitations with a reasonable expectation of success, in order to leverage the external device's existing communication hardware and processing capabilities to expand the digital hearing aid's functionality to audio signals generated or received by the external device itself, before or instead of signals received by the digital hearing aid itself. A POSITA would have been further motivated to modify this combined system, in view of Anderson or Jannard's disclosures, to include at least one rechargeable battery and/or configuring it to be able</p>

Claims	Prior Art
	<p>to be charged with and/or coupled to a docking station via a conductive element when placed within a charging apparatus to reduce power burden and size of the ear device, improve usability, and leverage known wireless and rechargeable power technologies known in the field at the time of the invention. A POSITA would have had a reasonable expectation of success since both the distributed processing architectures and rechargeable power systems were widely used and well-understood at the time, and their integration would have only required routine engineering to achieve predictable improvements in efficiency, down-time, and convenience.</p>
<p>[25b] wherein the rechargeable battery is able to be charged when the head-worn wearable electronic device is coupled to a docking station.</p>	<p>DHA systems, alone or in view of Anderson, expressly or inherently discloses “wherein the rechargeable battery is able to be charged when the head-worn wearable electronic device is coupled to a docking station.” Anderson ’783 at 22:11–61</p> <p>Stand-alone earpiece operation: If an earpiece cannot communicate with an RPU for any reason, stand-alone earpiece operation as a CIC hearing aid is generally desirable. Examples include occasions where the RPU must be turned off by pushing the RPU power-off key 78 (e.g., situations where operation of wireless transmitters is prohibited such as aircraft takeoffs and classified military briefings), RPU battery or other RPU electronics failure, misplacement of an RPU, and primary wireless link 17 interference. Note that a microphone located in the earpiece 10 is preferred to a body-worn microphone, since only an earpiece microphone 12 can be used to provide the desired default CIC hearing aid function (earpiece microphones also provide superior audio localization capabilities as the microphone moves with the user's head, and provide reduced wind and clothing noise compared with body-worn microphones). Parameters that control the stand-alone earpiece operation are stored in the earpiece CIC default parameter storage 870, which is accessed by the ASIC 845 to implement the desired stand-alone functions as described earlier. CIC default parameters include, but are not limited to, settings for the nominal and maximum allowable earpiece volume. Other information contained in the default parameter storage 870 includes, but is not limited to, the communication address and L/R (left/right) selection bit. To ensure that the CIC default parameter storage information cannot be changed inadvertently, changes are only allowed when the earpiece is not in use as, for example, when the earpiece battery 860 is being charged. A battery charging detector 840, that is sensitive to the high-level magnetic induction fields commonly used for earpiece battery charging, indicates to the ASIC 845 that the default parameter storage 870 values may be changed upon receipt of an appropriate command from the RPU. Specifically, to change information in the earpiece CIC default parameter storage 870, a user first activates the earpiece battery charging system, e.g., by placing the earpiece in a recharger/storage tray well known in the art. The RPU keyboard 946 is then used to enter a command to the DSP 948 that indicates the user's desire to change information in the default parameter storage 870, followed by the new information to be stored. The DSP 948 causes the RPU to transmit an interrogation with a reserved address code to which all earpieces being charged within</p>

Claims	Prior Art
	<p>the communication range of the RPU are receptive. The information that is to ultimately be loaded into the earpiece default parameter storage 870 is then transmitted by the RPU using the auxiliary bit field of many successive interrogations, said interrogations are received by the earpiece, and the desired information is loaded into the CIC default parameter storage 870 under control of the ASIC 845.</p> <p>A POSITA would have been motivated to combine the DHA system with the features of secondary references, including Anderson, to meet the claim limitations with a reasonable expectation of success for at least the reasons explained in claim 25.</p> <p>To the extent that the DHA system, alone or in view of Anderson, does not disclose any elements of this claim limitation, it would have been obvious to a person having ordinary skill in the art to combine, with the teachings of DHA systems, those elements, as known to them or as disclosed by other prior art identified in the cover pleading.</p>
<p>[26a] A personal audio system as recited in claim 24, wherein the head-worn wearable electronic device includes at least a rechargeable battery and at least one conductive element,</p>	<p>DHA systems, alone or in view of Anderson and Jannard, expressly or inherently discloses “[a] personal audio system as recited in claim 24, wherein the head-worn wearable electronic device includes at least a rechargeable battery and at least one conductive element.”</p> <p>Anderson ’783 at FIG. 8</p>  <p>The diagram, labeled FIG. 8, illustrates a personal audio system. At the top, an Antenna (40) is connected to an RF DEMOD & Modulator (805). The RF DEMOD & Modulator (805) is connected to a Video Amp with AGC (810), which is connected to a Video Pulse Quantizer (815), which in turn is connected to an Interrogation Data Register (820). A 60 MHz Oscillator (835) is connected to a mixer (830) and the RF DEMOD & Modulator (805). The mixer (830) is also connected to a Reply Data Register (825). The Interrogation Data Register (820) is connected to the Reply Data Register (825). The Reply Data Register (825) is connected to an ASIC (845). The ASIC (845) is connected to a Battery Charging Detector (840) and a Battery (860). The Battery (860) is connected to a 4-Bit A/D converter (865), which is connected to the ASIC (845). The ASIC (845) is connected to a Default Parameter Storage (870). The ASIC (845) is connected to a Microphone (875), which is connected to an 8-Bit Compressing A/D converter (880). The 8-Bit Compressing A/D converter (880) is connected to an 8-Bit Expanding D/A converter (850), which is connected to a Speaker (855).</p>

Claims	Prior Art
	<p>Anderson '783 at 3:40–41 FIG. 8 is a block diagram showing details of an earpiece for a full-featured preferred embodiment of the invention.</p> <p>Anderson '783 at 6:55–57 Earpiece power consumption is preferably minimized to extend battery life, and use of a rechargeable battery is desirable.</p> <p>Jannard '526 at FIG. 11</p>  <p style="text-align: center;">FIG. 11</p> <p>Jannard '526 at 4:27–28 FIG. 11 is a front elevational view of the audio device illustrated in FIG. 10.</p> <p>Jannard '526 at 24:56–25:3 The audio device 10D can also include power contacts 164 for recharging any rechargeable batteries connected thereto. For example, the power contacts 164 can be disposed on a lower edge of the orbitals 48D, 50D. Thus, with an appropriate recharging cradle (not shown), the audio device 10D can be laid on the cradle, thereby making contact between the power contacts 164 and corresponding contacts in the</p>

Claims	Prior Art
	<p>cradle (not shown). Alternatively, power contacts can be provided in numerous other locations as desired. For example, the power contacts 164 can be disposed at the ends of the ear stems 54D, 56D. A corresponding cradle can include two vertically oriented holes into which the ear stems are inserted for recharging. In this configuration, the lens within the orbitals 48D, 50D would face directly upwardly.</p> <p>Moreover, since the power contacts 164 make electrical contact with corresponding contacts in a recharging cradle to enable charging, the POSITA would have understood the power contacts 164 to be conductive elements.</p> <p>A POSITA would have been motivated to combine the teachings of the DHA system with Anderson and Jannard to meet the claim limitations with a reasonable expectation of success for the reasons described above with respect to claim [1a]. In addition, a POSITA would look to Jannard to supply these missing details—all of which were conventional—since both Anderson and Jannard disclose head-worn audio devices with rechargeable batteries, but Anderson does not expressly disclose the details of such a recharger.</p> <p>To the extent that the DHA system, alone or in view of Anderson and Abreu, does not disclose any elements of this claim limitation, it would have been obvious to a person having ordinary skill in the art to combine, with the teachings of DHA systems, those elements, as known to them or as disclosed by other prior art identified in the cover pleading.</p>
<p>[26b] wherein the rechargeable battery is able to be electrically coupled via the at least one conductive element to a charging apparatus and charged when the head-worn wearable electronic device is placed within the charging apparatus, and</p>	<p>DHA systems, alone or in view of Anderson and Jannard, expressly or inherently discloses “wherein the rechargeable battery is able to be electrically coupled via the at least one conductive element to a charging apparatus and charged when the head-worn wearable electronic device is placed within the charging apparatus.”</p> <p>To the extent that Anderson alone does not teach that a rechargeable battery is able to be electrically coupled via the at least one conductive element to the charging apparatus and charged when the head-worn wearable electronic device is placed within the charging apparatus, this element is met by Anderson in combination with Jannard.</p> <p>Jannard ’526 at 24:56–25:3</p> <p>The audio device 10D can also include power contacts 164 for recharging any rechargeable batteries connected thereto. For example, the power contacts 164 can be disposed on a lower edge of the orbitals 48D, 50D. Thus, with an appropriate recharging cradle (not shown), the audio device 10D can be laid on the cradle, thereby making contact between the power contacts 164 and corresponding contacts in the cradle (not shown). Alternatively, power contacts can be provided in numerous other locations as desired. For example, the power contacts 164 can be disposed at the ends of the ear stems 54D, 56D. A</p>

Claims	Prior Art
	<p>corresponding cradle can include two vertically oriented holes into which the ear stems are inserted for recharging. In this configuration, the lens within the orbitals 48D, 50D would face directly upwardly.</p> <p>A POSITA would have been motivated to combine the teachings of the DHA system with Anderson and Jannard to meet the claim limitations with a reasonable expectation of success for the reasons described above with respect to claim [1a] and [26a].</p> <p>To the extent that the DHA system, alone or in view of Anderson and Abreu, does not disclose any elements of this claim limitation, it would have been obvious to a person having ordinary skill in the art to combine, with the teachings of DHA systems, those elements, as known to them or as disclosed by other prior art identified in the cover pleading.</p>
<p>[26c] wherein the electrical coupling of the at least one conductive element with the charging apparatus is supported by at least a magnetic element.</p>	<p>DHA systems, alone or in view of Anderson and Jannard, expressly or inherently discloses “wherein the electrical coupling of the at least one conductive element with the charging apparatus is supported by at least a magnetic element.”</p> <p>For example, Jannard discloses that the power contacts 164 engage corresponding contacts in an “appropriate recharging cradle.” <i>Id.</i> at 24:59-63. To the extent not expressly disclosed in Jannard, it would have been obvious to the POSITA that the recharging cradle includes a magnetic element. This is because magnetic charging was well-known in the art by the effective filing date of the ’693 patent. It would have been obvious to use at least a magnetic element in Jannard’s charger to support electrical coupling of the power contacts to the corresponding contacts in the cradle in order to eliminate the need for charging cables.</p> <p>A POSITA would have been motivated to combine the teachings of the DHA system with Anderson and Jannard to meet the claim limitations with a reasonable expectation of success for the reasons described above with respect to claim [1a] and [26a].</p> <p>To the extent that the DHA system, alone or in view of Anderson and Abreu, does not disclose any elements of this claim limitation, it would have been obvious to a person having ordinary skill in the art to combine, with the teachings of DHA systems, those elements, as known to them or as disclosed by other prior art identified in the cover pleading.</p>
<p>27. A personal audio system as recited in claim 24, wherein the mobile communication device is configured to wirelessly communicate with the head-worn</p>	<p>DHA systems, alone or in view of Anderson, expressly or inherently discloses “[a] personal audio system as recited in claim 24, wherein the mobile communication device is configured to wirelessly communicate with the head-worn wearable electronic device over a wireless network.”</p>

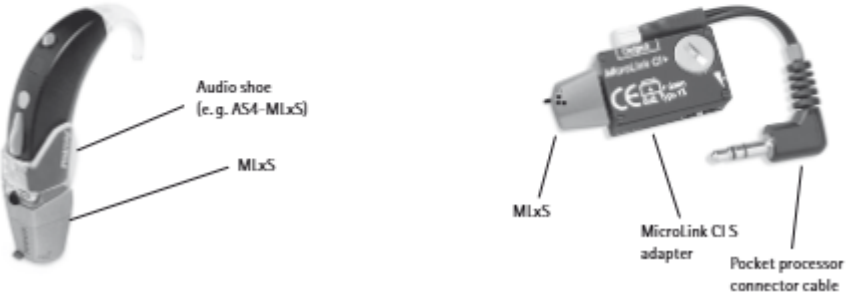
Claims	Prior Art
wearable electronic device over a wireless network.	<p data-bbox="583 175 1270 212">Phonak - SmartLink SX - FM Systems Manual at p. 3</p> <p data-bbox="638 228 1094 280">SmartLink SX is the first communication device to incorporate</p> <ul data-bbox="638 315 1031 521" style="list-style-type: none"> • Advanced FM features • Hearing Instrument remote control • A Bluetooth link for cell phone use <p data-bbox="583 548 1270 586">Phonak - SmartLink SX - FM Systems Manual at p. 4</p> <div data-bbox="638 602 1094 792"> <p data-bbox="638 602 932 678">Smart connections to all audio sources</p> <p data-bbox="638 683 1094 792">Any external audio input or an external microphone signal can be transmitted wirelessly via SmartLink to the hearing instrument, simply by attaching to the smart interface.</p> </div>

Claims	Prior Art
	<p data-bbox="583 175 1289 207">Phonak - SmartLink SX - FM Systems Manual at p. 10</p>  <p data-bbox="630 795 651 820">0</p> <p data-bbox="583 841 1289 873">Phonak - SmartLink SX - FM Systems Manual at p. 11</p> <p data-bbox="630 893 982 971">Smart digital synthesizer technology</p> <p data-bbox="630 974 1096 1334">Digital frequency synthesizer technology, only available with Phonak FM communication systems, enables online frequency synchronization and scanning of the frequencies used. Once a frequency is selected on the SmartLink, it automatically sends a signal to the FM receiver that adjusts the frequency of the receiver accordingly. This automatic synchronization process occurs every time the SmartLink is switched on. In cases where more than one FM transmitter is in use, at lectures or presentations for example, the SmartLink has an exclusive scan function. This seeks out the frequency being used by the other transmitter and sets the receiver frequency accordingly.</p>

Claims	Prior Art
	<p data-bbox="583 175 1289 212">Phonak - SmartLink SX - FM Systems Manual at p. 14</p>  <p data-bbox="583 841 1150 878">Phonak - SmartLink SX User Guide at p. 32</p> 

Claims	Prior Art
	<p>Phonak – Microlink User Guide at Page 2</p> <p>The desired signal – for example, a speaker’s voice or the sound from a TV – is picked up by the microphone which is placed near to the source of the sound. The transmitter sends an FM signal to the receiver, which immediately relays this to the hearing instrument for you to hear – crystal clear!</p> <p>Phonak – Annual Report (2003–2004) at p. 12</p> <p>Within the FM business segment, the Group continued to leverage existing key technology into new applications. With SmartLink, Phonak has achieved another industry first by allowing hearing instrument users to access other communications systems. SmartLink combines a remote control, a transmitter, a directional microphone system and a Bluetooth link to mobile phone all in one single device. Phonak also presented EduLink, an innovative FM device to assist the learning development of children with Auditory Processing Disorder (APD). EduLink was perceived very positively within the medical community and recent classroom tests show favourable results.</p> <p>Phonak – Online Interview re Perseo (April 7, 2003) at Page 9</p> <p><u>Liz Brassine/Phonak</u>: It is a digital microphone on the CIC. There are several advantages of a digital microphone over analog mics. First, increased dynamic range is provided. Additionally, with a digital mic, the analog-to-digital converter is housed within the microphone itself. This eliminates the wire between the AD-converter and the microphone. The wire is the source for ultra-high frequency interference, such as Bluetooth.</p> <p>Phonak – Claro 211 Product Information at Page 1</p> <p>FM features</p> <ul style="list-style-type: none"> – World’s smallest FM microchip – Compatible with all Phonak FM transmitters – Cordless TV and Radio link up – FM performance corresponding to audiological fitting – Soft Squelch system – Built in antenna – Band H (173–185 MHz); band J (186–201 MHz); band N (207–217 MHz)

Claims	Prior Art																		
	<div>Phonak – Claro 211 Product Information at Page 1</div> <div><div>Processing features</div><div><ul style="list-style-type: none">- Digital Perception Processing (DPP) in 20 critical bands- Gain and output controlled by DPP- 3 hearing programs (+ mute mode):<ul style="list-style-type: none">• QuietAdapt for hearing in relative quiet• NoiseAdapt for hearing speech in noise• Program 3 telephone default (other choices available)- FM programs: FM and FM+M- AutoSelect: automatic switching between QuietAdapt and NoiseAdapt- Claro's noise suppression technologies<ul style="list-style-type: none">• Adaptive digital AudioZoom (dAZ)• Fine-scale Noise Canceler in 20 bands• Soft squelch</div></div> <div>Phonak – Claro 211 Product Information at Page 2</div> <div><div>Principle of FM operation</div><div><p>The transmitter unit is worn by the person speaking, placed in front of the sound source or even connected directly to TV or HiFi. The transmitter unit collects the desired signal and transmits it via radio waves directly to the FM module integrated in the hearing instrument. The incoming FM signal received by the hearing instrument is processed either in accordance with the audiological setting of the instrument, providing consistent sound quality or according to an individual FM program.</p></div></div> <div>Phonak SmartLink SX Product Information at p. 2</div> <div><table><tr><th>Standards</th><th>Europe</th><th>USA</th></tr><tr><td>Wireless microphones</td><td>EN 300.422</td><td>FCC Part 15 Subpart C</td></tr><tr><td>Eletromagn. compatibility</td><td>EN 301.489</td><td>FCC Part 15 Subpart C</td></tr><tr><td>Short-range wireless transmission systems</td><td>EN 300.328</td><td>FCC Part 15 Subpart C</td></tr><tr><td>Safety regulations</td><td>EN 60950</td><td>FCC Part 15 Subpart C</td></tr><tr><td>Bluetooth</td><td colspan="2">*Bluetooth Special Interest Group* conformity test</td></tr></table></div>	Standards	Europe	USA	Wireless microphones	EN 300.422	FCC Part 15 Subpart C	Eletromagn. compatibility	EN 301.489	FCC Part 15 Subpart C	Short-range wireless transmission systems	EN 300.328	FCC Part 15 Subpart C	Safety regulations	EN 60950	FCC Part 15 Subpart C	Bluetooth	*Bluetooth Special Interest Group* conformity test	
Standards	Europe	USA																	
Wireless microphones	EN 300.422	FCC Part 15 Subpart C																	
Eletromagn. compatibility	EN 301.489	FCC Part 15 Subpart C																	
Short-range wireless transmission systems	EN 300.328	FCC Part 15 Subpart C																	
Safety regulations	EN 60950	FCC Part 15 Subpart C																	
Bluetooth	*Bluetooth Special Interest Group* conformity test																		

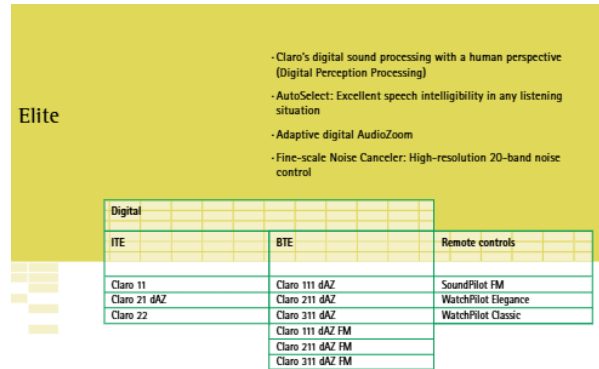
Claims	Prior Art
	<p>Phonak MicroLink MLxS Product Information at Page 1</p> <p>The MicroLink MLxS is the first miniaturized FM receiver with a digital frequency synthesizer for full multi-frequency flexibility. The receiver features fully automatic, non-contact frequency synchronization with either a WallPilot or the Campus S.</p> <p>The MicroLink MLxS is compatible with all Phonak FM transmitters and the WallPilot, and is recommended for use with this equipment.</p> <div data-bbox="638 441 1478 786"> <div data-bbox="638 441 945 461">MLxS with BTE instrument</div> <div data-bbox="1142 441 1457 461">MLxS with adapter for CI pocket processor</div>  </div> <p>Phonak – Media Release re Product Offerings in F2003/2004 at Page 3</p> <p>Phonak Communications introduces the new MicroLink receivers with frequency synthesizer</p> <p>Phonak Communications will introduce its new multi frequency MLx S and ML8 S receivers in the next few months. MLx S and ML8 S are the counterpart of Campus S transmitters, already introduced early 2002. Together, they comprise the new synthesized FM system generation, a system that transmits, without wires, speech to the hearing instrument in difficult listening situations.</p> <p>The key technology of the new system is its easy frequency adaptation. Given the extreme design constraints (low voltage, low power, miniature dimensions), the development of this technology represents an unprecedented breakthrough in wireless communication technology.</p>

Claims	Prior Art
	<p data-bbox="583 175 1409 207">Phonak and Other DHA – Manufacturers News Article at Page 1</p> <p data-bbox="636 224 945 516"><i>STÄFA, SWITZERLAND</i>—Phonak has introduced EduLink, which it predicts will contribute greatly to the successful treatment of children with auditory processing disorders (APDs). This new device is an "ear set" that transmits the teacher's voice, using FM technology, directly to the child's ear, thus effectively reducing ambient and background noise for improved speech discrimination. EduLink is designed to complement existing therapies in the treatment of APD.</p> <p data-bbox="583 540 1409 573">Phonak and Other DHA – Manufacturers News Article at Page 2</p> <p data-bbox="636 589 945 881">Phonak applies frequency synthesizer Phonak has integrated a new digital frequency synthesizer into several products: the Campus S transmitter, designed for use in learning environments; the MicroLink MLxS, the Perseo, the Claro dAZ FM S, and the Microlink CI S receiver for cochlear implant users. The digital synthesizer simplifies frequency management in schools and is intended to close the gap between children with and without hearing loss.</p> <p data-bbox="636 889 945 1076">This multi-frequency technology enables children with hearing instruments to move freely between classrooms, without the need for frequency re-setting. For wearers and teachers, the new FM technology means increased flexibility, ensuring that children and not technical issues are the focus of the education process.</p>

Claims

Prior Art

Phonak – Collections 2002 Product List at Page 5



Elite Line features:

- Clear, natural sound with precise, fully automatic volume control
- High level of user convenience with automatic hearing program selection
- FM receiver fully integrated in the BTE case (Claro Daz FM)
- Significantly improved speech intelligibility even in difficult listening situations
- Innovative high-tech solutions for Claro system control
- Design innovations for all system components
- State-of-the-art hearing computers for the optimised solution to a wide range of hearing problems

Claims

Prior Art

Phonak – Collections 2002 Product List at Page 9

BTE instruments	Claro 111 dAZ Claro 111 dAZ FM	Claro 211 dAZ Claro 211 dAZ FM	Claro 311 dAZ Claro 311 dAZ FM
			
			
Peak values G max./OSPL 90	52 dB/126 dB SPL	62 dB/135 dB SPL	77 dB/135 dB SPL
Features	Specially designed for steep high-frequency and mild hearing losses	Digital speech transparency	Digital speech transparency and higher gain
Description	Digital Perception Processing (DPP) AutoSelect (AS) Fine scale Noise Canceller (FNC) Adaptive digital AudioZoom (dAZ) FM module (111 dAZ FM) 5 hearing programs + "mute" Telecoil Battery size 13 Standard: TacTronic program selection switch Wind and weather protection	Digital Perception Processing (DPP) AutoSelect (AS) Fine scale Noise Canceller (FNC) Adaptive digital AudioZoom (dAZ) FM module (211 dAZ FM) 5 hearing programs + "mute" Telecoil Battery size 13 Standard: TacTronic program selection switch Wind and weather protection	Digital Perception Processing (DPP) AutoSelect (AS) Fine scale Noise Canceller (FNC) Adaptive digital AudioZoom (dAZ) FM module (311 dAZ FM) 5 hearing programs + "mute" Telecoil Battery size 13 Standard: TacTronic program selection switch Wind and weather protection
Options	Remote controls WatchPlot and SoundPlot FM Hooks: HE4, HE4 330, HE4 1000 Tamperproof battery cover ColorLine colors plus F25 (silver), F13 (milky transparent)	Remote controls WatchPlot and SoundPlot FM Hooks: HE4, HE4 330, HE4 1000 Tamperproof battery cover ColorLine colors plus F25 (silver), F13 (milky transparent)	Remote controls WatchPlot and SoundPlot FM Hooks: HE4, HE4 330, HE4 1000 Tamperproof battery cover ColorLine colors plus F25 (silver), F13 (milky transparent)
Programming	Loudness Perception Profile (LPP) with Active Response Unit (ARU) Fitting formulas: Phonak Digital, Phonak Digital Ski Slope and DSL (i/o)	Loudness Perception Profile (LPP) with Active Response Unit (ARU) Fitting formulas: Phonak Digital and DSL (i/o)	Loudness Perception Profile (LPP) with Active Response Unit (ARU) Fitting formulas: Phonak Digital and DSL (i/o)

Claims

Prior Art

Phonak – Collections 2002 Product List at Page 16

BTE instruments Digital technology	Acro 211*	Acro 211 AZ*
		
		
Peak values G max./OSPL 90	67 dB/129 dB SPL	67 dB/129 dB SPL
Features	Phonak's smallest digital BTE hearing computer	Phonak's smallest digital BTE hearing computer with Audio/zoom
Description	Digital MultiMode Signal Processing (dMSP) 3 hearing programs + "mute" High-resolution digital Noise Canceller (dNC) Feedback manager Customized signal tones for program confirmation and battery warning TacTronic multi-function switch Telecoil Audio-input/MicroLink compatible Battery size 13	Digital MultiMode Signal Processing (dMSP) 3 hearing programs + "mute" High-resolution digital Noise Canceller (dNC) Feedback manager Customized signal tones for program confirmation and battery warning TacTronic multi-function switch Telecoil Audio-input/MicroLink compatible Audio/zoom Battery size 13
Options	MicroLink ML5/MLx Audio shoe AS5 Tamperproof battery cover K55 ColorLine colors Hook: HES 1000 Mini ear hook	MicroLink ML5/MLx Audio shoe AS5 Tamperproof battery cover K55 ColorLine colors Hook: HES 1000 Mini ear hook

Claims

Prior Art

Phonak – Collections 2002 Product List at Page 17

Aero 311*/**	Aero 311 AZ*/**	Supero 412*	Supero 413 AZ*/**
			
  	  	  	  
77 dB/137 dB SPL	77 dB/137 dB SPL	84 dB/144 dB SPL	79 dB/144 dB SPL
High gain in a smaller case, digital BTE	High gain plus Audio/oom in a smaller case, digital BTE	Digital high power BTE in the high-end performance class	Digital high power BTE with unique Audio/oom multi-microphone technology
Digital MultiMode Signal Processing (dMSP) 3 hearing programs + "mute" High-resolution digital Noise Canceller (dNC) Feedback manager TacTronic multi-function switch Customised signal tones for program confirmation and battery warning Telecoil Audio-input/MicroLink compatible Battery size 13	Digital MultiMode Signal Processing (dMSP) 3 hearing programs + "mute" High-resolution digital Noise Canceller (dNC) Feedback manager TacTronic multi-function switch Customised signal tones for program confirmation and battery warning Telecoil Audio-input/MicroLink compatible Audio/oom Battery size 13	Digital MultiMode Signal Processing (dMSP) 5 hearing programs + "mute" Digital Noise Canceller (dNC) Power Dynamic Control for maximum protection without gain loss Feedback manager TacTronic multi-function switch Customised signal tones for program confirmation and battery warning Telecoil Audio-input/MicroLink compatible Battery size 6/5	Digital MultiMode Signal Processing (dMSP) 5 hearing programs + "mute" Digital Noise Canceller (dNC) Power Dynamic Control for maximum protection without gain loss Feedback manager TacTronic multi-function switch Customised signal tones for program confirmation and battery warning Telecoil Audio-input/MicroLink compatible Audio/oom Battery size 6/5
MicroLink M1.5/M1.x Audio shoe AS5 Tamperproof battery cover KS5 ColorLine colors Hooks: HE2, HE2 680, HE2 1000 Mini ear hooks	MicroLink M1.5/M1.x Audio shoe AS5 Tamperproof battery cover KS5 ColorLine colors Hooks: HE2, HE2 680, HE2 1000 Mini ear hooks	MicroLink M1.4/M1.x Audio shoe AS4 Tamperproof battery cover KS 5-3 ColorLine colors Hooks: HE3, HE3 680, HE3 1000 Remote control SoundPilot FM Remote control WatchPilot	MicroLink M1.4/M1.x Audio shoe AS4 Tamperproof battery cover KS 5-3 ColorLine colors Hooks: HE3, HE3 680, HE3 1000 Remote control SoundPilot FM Remote control WatchPilot

Claims

Prior Art

Phonak – Collections 2002 Product List at Page 18

Digitally programmable BTE instruments	Piconet2 P2	Piconet2 P2 AZ
		
		
Peak values G max./OSPL 90	62 dB/128 dB SPL	62 dB/128 dB SPL
Features	Phonak's smallest digitally programmable BTE hearing computer	Phonak's smallest digitally programmable BTE hearing computer with Audio/oom
Description	Multi Dynamic Compression Control (MDCC) 3 hearing programs Telecoil Audio-input/MicroLink compatible Battery size 13	Multi Dynamic Compression Control (MDCC) 3 hearing programs Telecoil Audio-input/MicroLink compatible Audio/oom (also available without DHC) Battery size 13
Options	MicroLink ML6/MLx Audio shoe AS6 Tamperproof battery cover KS5 ColorLine colors Hooks: HE2 1000, HE5 1000 Mini ear hooks Remote control: DHC-2, DHC-4	MicroLink ML6/MLx Audio shoe AS6 Tamperproof battery cover KS5 ColorLine colors Hooks: HE2 1000, HE5 1000 Mini ear hooks Remote control: DHC-2, DHC-4

Claims	Prior Art																												
	<p>Phonak – Collections 2002 Product List at Page 19</p> <table><thead><tr><th>SonoForte2 P3 AZ</th><th>PowerZoom P4 AZ</th><th>NovoForte E3</th><th>NovoForte E4</th></tr></thead><tbody><tr><td></td><td></td><td></td><td></td></tr><tr><td></td><td></td><td></td><td></td></tr><tr><td>77 dB/137 dB SPL</td><td>79 dB/144 dB SPL</td><td>77 dB/137 dB SPL</td><td>84 dB/144 dB SPL</td></tr><tr><td>High gain plus Audio/oom in a smaller case</td><td>High gain and Audio/oom</td><td>2 switch-selectable programs</td><td>Phonak's most powerful digitally programmable BTE instrument</td></tr><tr><td>Multi Dynamic Compression Control (MDCC) 3 hearing programs Telecoil Audio-input/Microlink compatible Audio/oom (also available without DHC) Battery size 13</td><td>Multi Dynamic Compression Control (MDCC) 3 hearing programs Telecoil Audio-input/Microlink compatible Audio/oom (also available without DHC) Battery size 6/5</td><td>Signal processing strategy: WDRC, SC + a.R.T., PC 2 switch-selectable hearing programs Telecoil Audio-input/Microlink compatible Programmable user controls Battery size 6/5</td><td>Signal processing strategy: WDRC, SC + a.R.T., PC 2 switch-selectable hearing programs Telecoil Audio-input/Microlink compatible Programmable user controls Battery size 6/5</td></tr><tr><td>Microlink ML6/MLx Audio shoe AS6 Tamperproof battery cover KS5 Colorline colors Hooks: HE2, HE2 680, HE2 1000 Mini ear hooks Remote control: DHC-2, DHC-4</td><td>Microlink ML7/MLx Audio shoe AS7 Tamperproof battery cover KS7 Colorline colors Hooks: HE3, HE3 680, HE3 1000 Mini ear hooks Remote control: DHC-2, DHC-4</td><td>Microlink ML7/MLx Audio shoe AS7 Tamperproof battery cover KS7 Colorline colors Hooks: HE3, HE3 680, HE3 1000 Mini ear hooks Remote control: DHC-3</td><td>Microlink ML7/MLx Audio shoe AS7 Tamperproof battery cover KS7 Colorline colors Hooks: HE3, HE3 680, HE3 1000 Mini ear hooks Remote control: DHC-3</td></tr></tbody></table>	SonoForte2 P3 AZ	PowerZoom P4 AZ	NovoForte E3	NovoForte E4									77 dB/137 dB SPL	79 dB/144 dB SPL	77 dB/137 dB SPL	84 dB/144 dB SPL	High gain plus Audio/oom in a smaller case	High gain and Audio/oom	2 switch-selectable programs	Phonak's most powerful digitally programmable BTE instrument	Multi Dynamic Compression Control (MDCC) 3 hearing programs Telecoil Audio-input/Microlink compatible Audio/oom (also available without DHC) Battery size 13	Multi Dynamic Compression Control (MDCC) 3 hearing programs Telecoil Audio-input/Microlink compatible Audio/oom (also available without DHC) Battery size 6/5	Signal processing strategy: WDRC, SC + a.R.T., PC 2 switch-selectable hearing programs Telecoil Audio-input/Microlink compatible Programmable user controls Battery size 6/5	Signal processing strategy: WDRC, SC + a.R.T., PC 2 switch-selectable hearing programs Telecoil Audio-input/Microlink compatible Programmable user controls Battery size 6/5	Microlink ML6/MLx Audio shoe AS6 Tamperproof battery cover KS5 Colorline colors Hooks: HE2, HE2 680, HE2 1000 Mini ear hooks Remote control: DHC-2, DHC-4	Microlink ML7/MLx Audio shoe AS7 Tamperproof battery cover KS7 Colorline colors Hooks: HE3, HE3 680, HE3 1000 Mini ear hooks Remote control: DHC-2, DHC-4	Microlink ML7/MLx Audio shoe AS7 Tamperproof battery cover KS7 Colorline colors Hooks: HE3, HE3 680, HE3 1000 Mini ear hooks Remote control: DHC-3	Microlink ML7/MLx Audio shoe AS7 Tamperproof battery cover KS7 Colorline colors Hooks: HE3, HE3 680, HE3 1000 Mini ear hooks Remote control: DHC-3
SonoForte2 P3 AZ	PowerZoom P4 AZ	NovoForte E3	NovoForte E4																										
																													
																													
77 dB/137 dB SPL	79 dB/144 dB SPL	77 dB/137 dB SPL	84 dB/144 dB SPL																										
High gain plus Audio/oom in a smaller case	High gain and Audio/oom	2 switch-selectable programs	Phonak's most powerful digitally programmable BTE instrument																										
Multi Dynamic Compression Control (MDCC) 3 hearing programs Telecoil Audio-input/Microlink compatible Audio/oom (also available without DHC) Battery size 13	Multi Dynamic Compression Control (MDCC) 3 hearing programs Telecoil Audio-input/Microlink compatible Audio/oom (also available without DHC) Battery size 6/5	Signal processing strategy: WDRC, SC + a.R.T., PC 2 switch-selectable hearing programs Telecoil Audio-input/Microlink compatible Programmable user controls Battery size 6/5	Signal processing strategy: WDRC, SC + a.R.T., PC 2 switch-selectable hearing programs Telecoil Audio-input/Microlink compatible Programmable user controls Battery size 6/5																										
Microlink ML6/MLx Audio shoe AS6 Tamperproof battery cover KS5 Colorline colors Hooks: HE2, HE2 680, HE2 1000 Mini ear hooks Remote control: DHC-2, DHC-4	Microlink ML7/MLx Audio shoe AS7 Tamperproof battery cover KS7 Colorline colors Hooks: HE3, HE3 680, HE3 1000 Mini ear hooks Remote control: DHC-2, DHC-4	Microlink ML7/MLx Audio shoe AS7 Tamperproof battery cover KS7 Colorline colors Hooks: HE3, HE3 680, HE3 1000 Mini ear hooks Remote control: DHC-3	Microlink ML7/MLx Audio shoe AS7 Tamperproof battery cover KS7 Colorline colors Hooks: HE3, HE3 680, HE3 1000 Mini ear hooks Remote control: DHC-3																										

Claims	Prior Art																								
	<div>Phonak – Collections 2002 Product List at Page 24</div> <div><div><div>Supero 411*</div><div></div></div><div><div>PicoForte3 SC-D</div><div></div></div><div><div>PicoForte3 PP-SC</div><div></div></div></div> <div>BTE instruments</div> <table><tr><th></th><th></th><th></th><th></th></tr><tr><td>Peak values G max./OSPL 90</td><td>79 dB/142 dB SPL</td><td>64 dB/127 dB SPL</td><td>78 dB/133 dB SPL</td></tr><tr><td>Features</td><td>Digital high-power BTE</td><td>Directional microphone, switchable</td><td>High power</td></tr><tr><td>Description</td><td>d.limiting signal processing strategy 1 hearing program + "mute" Power Dynamic Control for maximum protection without gain loss Feedback manager TacTronic multi-function switch Customised signal tones for program confirmation and battery warning Telecoil Audio-input/Microlink compatible Battery size 675</td><td>Super Compression with a.R.L. Telecoil Audio-input/Microlink compatible Battery size 13</td><td>Super Compression + a.R.L. Telecoil Audio-input/Microlink compatible Battery size 13</td></tr><tr><td>Options</td><td>Microlink ML4/MLx Audio shoe ASA Tamperproof battery cover KS 5-3 Colorline colors Hooks: HE3, HE3 GR0, H3 1000</td><td>Microlink ML5/MLx Audio shoe ASSA Tamperproof battery cover KS 5-3 Colorline colors Hooks: HE2, HE2 GR0, H2 1000 Mini ear hooks</td><td>Microlink ML5/MLx Audio shoe ASSA Tamperproof battery cover KS 5-3 Colorline colors Hooks: HE2, HE2 GR0, H2 1000 Mini ear hooks</td></tr><tr><td>Programming</td><td>Programmable with PC (IBM compatible) and HI-PRO programming interface PFG Software, version B.0 or above</td><td>LC: Low cut, 6 levels SSPL: Output, 6 levels Fitting recommendation by PFG based on audiogram</td><td>LC: Low cut, 6 levels SSPL: Output, 6 levels Fitting recommendation by PFG based on audiogram</td></tr></table>					Peak values G max./OSPL 90	79 dB/142 dB SPL	64 dB/127 dB SPL	78 dB/133 dB SPL	Features	Digital high-power BTE	Directional microphone, switchable	High power	Description	d.limiting signal processing strategy 1 hearing program + "mute" Power Dynamic Control for maximum protection without gain loss Feedback manager TacTronic multi-function switch Customised signal tones for program confirmation and battery warning Telecoil Audio-input/Microlink compatible Battery size 675	Super Compression with a.R.L. Telecoil Audio-input/Microlink compatible Battery size 13	Super Compression + a.R.L. Telecoil Audio-input/Microlink compatible Battery size 13	Options	Microlink ML4/MLx Audio shoe ASA Tamperproof battery cover KS 5-3 Colorline colors Hooks: HE3, HE3 GR0, H3 1000	Microlink ML5/MLx Audio shoe ASSA Tamperproof battery cover KS 5-3 Colorline colors Hooks: HE2, HE2 GR0, H2 1000 Mini ear hooks	Microlink ML5/MLx Audio shoe ASSA Tamperproof battery cover KS 5-3 Colorline colors Hooks: HE2, HE2 GR0, H2 1000 Mini ear hooks	Programming	Programmable with PC (IBM compatible) and HI-PRO programming interface PFG Software, version B.0 or above	LC: Low cut, 6 levels SSPL: Output, 6 levels Fitting recommendation by PFG based on audiogram	LC: Low cut, 6 levels SSPL: Output, 6 levels Fitting recommendation by PFG based on audiogram
																									
Peak values G max./OSPL 90	79 dB/142 dB SPL	64 dB/127 dB SPL	78 dB/133 dB SPL																						
Features	Digital high-power BTE	Directional microphone, switchable	High power																						
Description	d.limiting signal processing strategy 1 hearing program + "mute" Power Dynamic Control for maximum protection without gain loss Feedback manager TacTronic multi-function switch Customised signal tones for program confirmation and battery warning Telecoil Audio-input/Microlink compatible Battery size 675	Super Compression with a.R.L. Telecoil Audio-input/Microlink compatible Battery size 13	Super Compression + a.R.L. Telecoil Audio-input/Microlink compatible Battery size 13																						
Options	Microlink ML4/MLx Audio shoe ASA Tamperproof battery cover KS 5-3 Colorline colors Hooks: HE3, HE3 GR0, H3 1000	Microlink ML5/MLx Audio shoe ASSA Tamperproof battery cover KS 5-3 Colorline colors Hooks: HE2, HE2 GR0, H2 1000 Mini ear hooks	Microlink ML5/MLx Audio shoe ASSA Tamperproof battery cover KS 5-3 Colorline colors Hooks: HE2, HE2 GR0, H2 1000 Mini ear hooks																						
Programming	Programmable with PC (IBM compatible) and HI-PRO programming interface PFG Software, version B.0 or above	LC: Low cut, 6 levels SSPL: Output, 6 levels Fitting recommendation by PFG based on audiogram	LC: Low cut, 6 levels SSPL: Output, 6 levels Fitting recommendation by PFG based on audiogram																						

Claims

Prior Art

Phonak – Collections 2002 Product List at Page 25

PicoFort3 PP-C-P



PicoFort3 PP-C-L-P



SuperFront PP-C-2



SuperFront PP-C-4



Other products 3
Superfront

FM	FM	FM	FM
84 dB/137 dB SPL	82 dB/137 dB SPL	79 dB/142 dB SPL	82 dB/142 dB SPL
High gain	Extra power in low-frequency range	High power instrument	High-power instrument
Push-pull amplifier Peak Clipping Telecoil Audio input/Microlink compatible Battery size 13	Push-pull amplifier Peak Clipping Telecoil Audio input/Microlink compatible Battery size 13	Push-pull amplifier Peak Clipping Telecoil Audio input/Microlink compatible Battery size 675	Push-pull amplifier Peak Clipping Telecoil Audio input/Microlink compatible Battery size 675
Microlink M5/Mix Audio shoe AS5A Tamperproof battery cover KS 5-3 ColorLine colors Hooks: H12, H12 G80, H2 1000 Mini ear hooks	Microlink M5/Mix Audio shoe AS5A Tamperproof battery cover KS 5-3 ColorLine colors Hooks: H12, H12 G80, H2 1000 Mini ear hooks	Microlink M5/Mix Audio shoe AS3A Tamperproof battery cover KS3 ColorLine colors Hooks: H12, H12 G80, PD 1000 Mini ear hooks	Microlink M5/Mix Audio shoe AS3A Tamperproof battery cover KS3 ColorLine colors Hooks: H12, H12 G80, PD 1000 Mini ear hooks
IC: Low cut, 6 levels SSPL: Output, 6 levels Fitting recommendation by PFG based on audiogram	IC: Low cut, 6 levels SSPL: Output, 6 levels Fitting recommendation by PFG based on audiogram	IC: Low cut, 5 levels SSPL: Output, 7 levels Fitting recommendation by PFG based on audiogram	IC: Low cut, 9 levels HQB: High cut, 3 levels; bypass, 4 levels SSPL: Output, 6 levels G: Gain, 6 levels Fitting recommendation by PFG based on audiogram

Claims

Prior Art

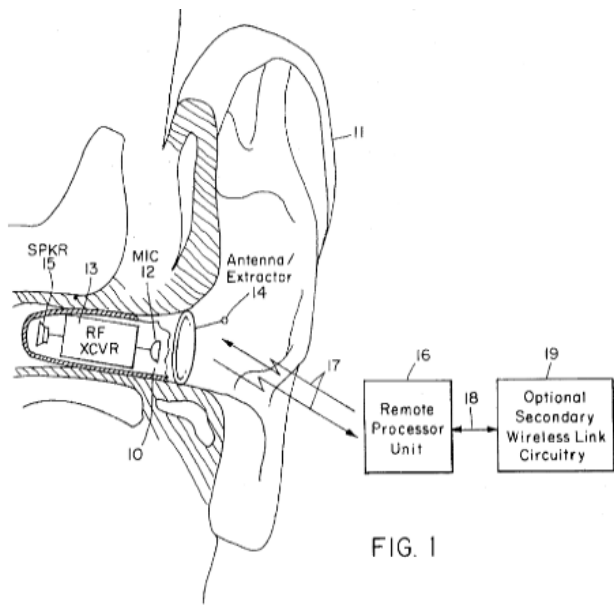
Phonak – Collections 2002 Product List at Page 26

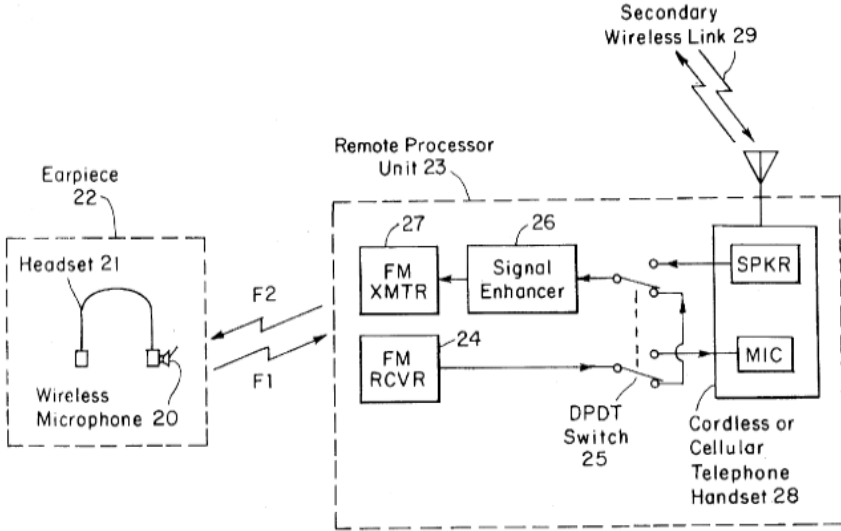
	SuperFront PP-C-L-4 SuperFront PP-C-L-4+	SuperFront PP-SC	Classica AGC
BTE instruments			
			
Peak values G max./OSPL 90	87 dB/142 dB SPL (PP-C-L-4) 81 dB/142 dB SPL (PP-C-L-4+)	87 dB/140 dB SPL	65 dB/133 dB SPL
Features	High-power low-freq. range (PP-C-L-4) Super low-freq. instrument (PP-C-L-4+)	High power	AGC: 0 with a.R.I.
Description	Push-pull amplifier Peak Clipping Telecoil Audio input/MicroLink compatible Battery size 675	Push-pull amplifier Super Compression with a.R.I. Telecoil Audio input/MicroLink compatible Battery size 675	Broadband response Telecoil Audio input/MicroLink compatible Battery size 675
Options	MicroLink M1.3/MLx Audio shoe AS2 Tamperproof battery cover KS3 Colorline colors Hooks: HE, HE 680, PD 1000 Mini ear hooks	MicroLink M1.3/MLx Audio shoe AS2 Tamperproof battery cover KS3 Colorline colors Hooks: HE, HE 680, PD 1000 Mini ear hooks	MicroLink M1.4/MLx Audio shoe AS4 Tamperproof battery cover KS4 Colorline colors Hooks: HE3, HE3 680, HE3 1000 Mini ear hooks
Programming	IC: Low cut, 9 levels HC: High cut, 6 levels SSPL: Output, 6 levels G: Gain, 6 levels Fitting recommendation by PFG based on audiogram	IC: Low cut, 9 levels BP/RC: Bypass, 6 levels; high cut, 3 levels SSPL: Output, 9 levels G: Gain, 6 levels Fitting recommendation by PFG based on audiogram	IC: Low cut, 5 levels SSPL: Output, 5 levels Fitting recommendation by PFG based on audiogram

Claims	Prior Art																																
	<p>Phonak – Collections 2002 Product List at Page 27</p> <table><thead><tr><th>Classica C-D</th><th>Classica PP-SC</th><th>Classica PP-C-P</th><th>Classica PP-C-L-P</th></tr></thead><tbody><tr><td></td><td></td><td></td><td></td></tr><tr><td> FM</td><td> FM</td><td> FM</td><td> FM</td></tr><tr><td>70 dB/139 dB SPL</td><td>80 dB/137 dB SPL</td><td>84 dB/144 dB SPL</td><td>80 dB/144 dB SPL</td></tr><tr><td>Broadband BTE with directional microphone</td><td>High power</td><td>High gain</td><td>Extra power in low frequency range</td></tr><tr><td>Push-pull amplifier Peak Clipping Telecoil Audio-input/Microlink compatible Battery size 6/75</td><td>Push-pull amplifier Super Compression with a.R.L. Telecoil Audio-input/Microlink compatible Battery size 6/75</td><td>Push-pull amplifier Peak Clipping Broadband Telecoil Audio-input/Microlink compatible Battery size 6/75</td><td>Push-pull amplifier Peak Clipping Broadband Telecoil Audio-input/Microlink compatible Battery size 6/75</td></tr><tr><td>Microlink M14/MLx Audio shoe AS4 Tamperproof battery cover KS4 ColorLine colors Hooks: HE3, HE3 680, HE3 1000 Mini ear hooks</td><td>Microlink M14/MLx Audio shoe AS4 Tamperproof battery cover KS4 ColorLine colors Hooks: HE3, HE3 680, HE3 1000 Mini ear hooks</td><td>Microlink M14/MLx Audio shoe AS4 Tamperproof battery cover KS4 ColorLine colors Hooks: HE3, HE3 680, HE3 1000 Mini ear hooks</td><td>Microlink M14/MLx Audio shoe AS4 Tamperproof battery cover KS4 ColorLine colors Hooks: HE3, HE3 680, HE3 1000 Mini ear hooks</td></tr><tr><td>IC: Low cut, 5 levels SSPL: Output, 5 levels Fitting recommendation by PFG based on audiogram</td><td>IC: Low cut, 5 levels SSPL: Output, 5 levels Fitting recommendation by PFG based on audiogram</td><td>IC: Low cut, 5 levels SSPL: Output, 5 levels Fitting recommendation by PFG based on audiogram</td><td>IC: Low cut, 5 levels SSPL: Output, 5 levels Fitting recommendation by PFG based on audiogram</td></tr></tbody></table>	Classica C-D	Classica PP-SC	Classica PP-C-P	Classica PP-C-L-P					 FM	 FM	 FM	 FM	70 dB/139 dB SPL	80 dB/137 dB SPL	84 dB/144 dB SPL	80 dB/144 dB SPL	Broadband BTE with directional microphone	High power	High gain	Extra power in low frequency range	Push-pull amplifier Peak Clipping Telecoil Audio-input/Microlink compatible Battery size 6/75	Push-pull amplifier Super Compression with a.R.L. Telecoil Audio-input/Microlink compatible Battery size 6/75	Push-pull amplifier Peak Clipping Broadband Telecoil Audio-input/Microlink compatible Battery size 6/75	Push-pull amplifier Peak Clipping Broadband Telecoil Audio-input/Microlink compatible Battery size 6/75	Microlink M14/MLx Audio shoe AS4 Tamperproof battery cover KS4 ColorLine colors Hooks: HE3, HE3 680, HE3 1000 Mini ear hooks	Microlink M14/MLx Audio shoe AS4 Tamperproof battery cover KS4 ColorLine colors Hooks: HE3, HE3 680, HE3 1000 Mini ear hooks	Microlink M14/MLx Audio shoe AS4 Tamperproof battery cover KS4 ColorLine colors Hooks: HE3, HE3 680, HE3 1000 Mini ear hooks	Microlink M14/MLx Audio shoe AS4 Tamperproof battery cover KS4 ColorLine colors Hooks: HE3, HE3 680, HE3 1000 Mini ear hooks	IC: Low cut, 5 levels SSPL: Output, 5 levels Fitting recommendation by PFG based on audiogram	IC: Low cut, 5 levels SSPL: Output, 5 levels Fitting recommendation by PFG based on audiogram	IC: Low cut, 5 levels SSPL: Output, 5 levels Fitting recommendation by PFG based on audiogram	IC: Low cut, 5 levels SSPL: Output, 5 levels Fitting recommendation by PFG based on audiogram
Classica C-D	Classica PP-SC	Classica PP-C-P	Classica PP-C-L-P																														
																																	
 FM	 FM	 FM	 FM																														
70 dB/139 dB SPL	80 dB/137 dB SPL	84 dB/144 dB SPL	80 dB/144 dB SPL																														
Broadband BTE with directional microphone	High power	High gain	Extra power in low frequency range																														
Push-pull amplifier Peak Clipping Telecoil Audio-input/Microlink compatible Battery size 6/75	Push-pull amplifier Super Compression with a.R.L. Telecoil Audio-input/Microlink compatible Battery size 6/75	Push-pull amplifier Peak Clipping Broadband Telecoil Audio-input/Microlink compatible Battery size 6/75	Push-pull amplifier Peak Clipping Broadband Telecoil Audio-input/Microlink compatible Battery size 6/75																														
Microlink M14/MLx Audio shoe AS4 Tamperproof battery cover KS4 ColorLine colors Hooks: HE3, HE3 680, HE3 1000 Mini ear hooks	Microlink M14/MLx Audio shoe AS4 Tamperproof battery cover KS4 ColorLine colors Hooks: HE3, HE3 680, HE3 1000 Mini ear hooks	Microlink M14/MLx Audio shoe AS4 Tamperproof battery cover KS4 ColorLine colors Hooks: HE3, HE3 680, HE3 1000 Mini ear hooks	Microlink M14/MLx Audio shoe AS4 Tamperproof battery cover KS4 ColorLine colors Hooks: HE3, HE3 680, HE3 1000 Mini ear hooks																														
IC: Low cut, 5 levels SSPL: Output, 5 levels Fitting recommendation by PFG based on audiogram	IC: Low cut, 5 levels SSPL: Output, 5 levels Fitting recommendation by PFG based on audiogram	IC: Low cut, 5 levels SSPL: Output, 5 levels Fitting recommendation by PFG based on audiogram	IC: Low cut, 5 levels SSPL: Output, 5 levels Fitting recommendation by PFG based on audiogram																														

Claims	Prior Art
	<p data-bbox="583 175 1234 207">Phonak – Collections 2002 Product List at Page 30</p> <div data-bbox="625 224 1323 824">  <p data-bbox="940 228 1318 321">MLx, the most flexible FM receiver ever made, can be used not only with Phonak BTEs but also with most third-party hearing instruments available on the market today. This is possible thanks to the standard Euro connector which simply plugs into the audio shoe of the hearing instrument.</p> <p data-bbox="940 342 1276 375">Once again, Phonak demonstrates its commitment to innovative developments that serve hearing-impaired adults and children.</p> <p data-bbox="632 342 724 375">MLx compatible BTE makes:</p> <ul data-bbox="632 412 808 820" style="list-style-type: none"> Argosy Audifon Belton Danavox GN ReSound Hansaton Interton Oticon Philips Phonak Rexton Siemens Sonic Innovations Unitron Widex <p data-bbox="842 639 1050 721">Detailed, up-to-date compatibility lists of Phonak's MLx MicroLink receiver are available online at: www.phonak.com/mlx</p> </div> <p data-bbox="583 846 1612 878">Lewis – Speech perception in noise: directional microphones (2004) at Page 429</p> <p data-bbox="632 894 1024 1019">FM technology was used. As previously mentioned, Phonak Corporation developed a new personal FM system receiver, the Phonak Microlink, which does not utilize wires or a body-worn box like its predecessors. This</p>

Claims	Prior Art
	<p data-bbox="583 175 1608 211">Lewis – Speech perception in noise: directional microphones (2004) at Page 431</p> <p data-bbox="636 228 1016 816"> All subjects were fit with digital Phonak Claro 311 dAZ BTE hearing aids bilaterally. All earmolds had select-a-vent (SAV) venting and #13 or 3 mm horn tubing. In addition to the hearing aids, subjects were fit with Phonak Microlink ML8 FM receivers bilaterally. These FM receivers attach to the bottom of a BTE hearing aid and may be utilized in either the “FM only” mode, which attenuates the hearing aid microphone by 20 dB, or in the FM plus hearing aid mode, which allows for FM input and input of environmental sounds via the hearing aid microphones simultaneously at the same output level. The Phonak TX3 HandyMic FM transmitter served as the FM transmitter. This transmitter has three microphone options: (1) Wide Angle, which picks up sounds arriving from all directions around the transmitter microphone equally; (2) Zoom, which provides reduced amplification to signals arriving from the rear (cardioid); and (3) SuperZoom, which provides reduced amplification to signals originating from the rear and the sides (hypercardioid). The </p> <p data-bbox="583 842 1232 878">Oticon – DigiFocus White Paper (1996) at Page 11</p> <p data-bbox="636 896 1308 1029"> Many people find it uncomfortable to be wired to a computer. This “constricted” feeling can be avoided by using the OtiSet programming belt, which is worn on the shoulder. The programming belt provides cordless infrared communication between the PC and the hearing instrument while in the ear. Alternatively, you may use the HI-PRO interface for a wired solution. </p> <p data-bbox="636 1052 1308 1209"> OtiSet offers a choice of three different views of the fitting. The Audiological View provides a visual guide to instrument performance in four typical sound environments. The Technical View illustrates traditional technical information such as coupler data and insertion gain at different input levels. With the Binaural View you are able to examine settings for both right and left instruments at the same time. </p> <p data-bbox="583 1235 924 1268">Anderson ’783 at Abstract</p> <p data-bbox="636 1274 1990 1450"> A hearing aid or audio communication system includes an earpiece (10) that can be hidden in the ear canal, and which communicates wirelessly with a remote processor unit, or RPU (16), that enhances audio signals and can be concealed under clothing. Sounds from the environment are picked up by a microphone (12) in the earpiece and sent with other information over a two-way wireless link (17) to the RPU (16). The wireless link (17) uses microwaves for component miniaturization. Furthermore, use of radar </p>

Claims	Prior Art
	<p>technology to implement the wireless link (17), with an RPU (16) interrogator and earpiece (10) transponder, reduces earpiece size and power, as no microwave oscillator is needed in the earpiece (10). Optional secondary wireless link circuitry (19) can be used between the RPU (16) and a cellular telephone system or other sources of information. Electronic voice recognition and response can control system operation.</p> <p>Anderson '783 at FIG. 1</p>  <p>FIG. 1</p>

Claims	Prior Art
	<p data-bbox="583 175 898 207">Anderson '783 at FIG. 2</p>  <p data-bbox="1024 789 1117 821">FIG. 2</p> <p data-bbox="583 846 905 878">Anderson '783 at 1:5–18</p> <p data-bbox="627 883 1953 1133">The present invention is directed to hearing aids, and in particular to a hearing aid having an earpiece housing worn in or at the ear and a remote processor unit (RPU) worn by or located near the user that wirelessly receives signals from and transmits signals to the earpiece. The present invention is also directed to the use of hearing aid systems as wireless communication devices for hands-free cellular telephone and mobile radio communication "handsets," covert operation and control of hearing aids, hearing protection and noise cancellation simultaneous with binaural hearing aid functions, hearing test equipment, location of misplaced hearing aid system components and wireless cochlear implants.</p> <p data-bbox="583 1190 919 1222">Anderson '783 at 1:50–65</p> <p data-bbox="627 1227 1984 1442">The present invention uses an audio path that is different from the audio path used by devices in the prior art. The present invention uses a primary two-way wireless link (not a one-way link as in the prior art) between an earpiece worn in or at the user's ear and an RPU worn by or located near the user. In the present invention, audio signals from the environment are picked up by a microphone in the earpiece (not a hand-held microphone as used by wireless hearing aids in the prior art) and transmitted over the primary two-way wireless link to the RPU (instead of going to processor electronics contained in the earpiece, as</p>

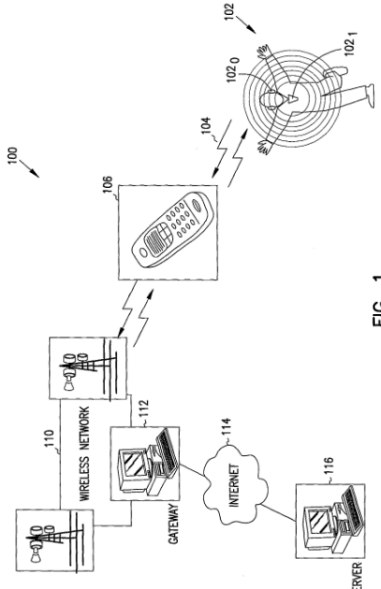
Claims	Prior Art
	<p>in basic hearing aids known in the prior art), where the audio signals are enhanced according to the user's needs before transmission over the primary wireless link to the earpiece. Signal processing is performed in the RPU rather than the earpiece to take advantage of relaxed size and power constraints.</p> <p>Anderson '783 at 2:20–39</p> <p>It is an object of the present invention to provide a new and useful auditory aid for hearing-impaired persons (i.e., those having certain residual hearing) by removing audio signal enhancement functions from the earpiece and placing them in an RPU. Use of an RPU provides several advantages over systems that attempt to place all system capabilities within the earpiece. The RPU approach allows a simple earpiece design comprising a miniature low-power wireless transceiver, microphone and speaker. Note that the speaker is also known in the hearing aid art as a "receiver," but the term "speaker" is used here to avoid confusion (similarly, the term "talker," not "speaker," is used to describe a human producing vocal sounds). The resulting earpiece is extremely small, can be hidden from view in the ear canal if desired, and allows complete freedom of movement when a primary two-way wireless link to the RPU is used. Processing for all major system capabilities, such as amplification and other forms of signal enhancement, takes place in the RPU where size and power constraints are relaxed, leading to a cost-effective design.</p> <p>Anderson '783 at 2:60–3:5</p> <p>A further object of this invention is to provide some degree of protection for the residual hearing capability of a hearing-impaired user (or a non-impaired user who wishes to avoid impairment) in a wireless hearing aid system having the simultaneous capabilities of noise cancellation and binaural processing (e.g., directionality). Other objects of the invention are to provide a convenient means for testing a user's hearing capability without the need for additional equipment, to assist users in the location of misplaced hearing aid system components, and to provide profoundly deaf cochlear implant patients with a wireless system allowing improved freedom of movement compared to existing wire-connected systems.</p> <p>Anderson '783 at 3:52–4:14</p> <p>One preferred embodiment of the invention herein described is shown in block diagram form in FIG. 1. In this figure, the earpiece 10 shown worn in the ear 11 uses a standard completely-in-the-canal (CIC) housing, well known in the art, although many other earpiece housing types known in the art (e.g., behind-the-ear, or BTE) can be used. The earpiece shown in FIG. 1 comprises a microphone 12, RF transceiver 13 with an antenna that doubles as an earpiece extractor 14, and speaker 15. Although the dual-use antenna/extractor arrangement shown in FIG. 1 is often desirable, a separate antenna or antennae</p>

Claims	Prior Art
	<p>can be used in conjunction with many different types of earpiece extractors as appropriate for a particular application. Antennae and extractors can be disguised as ear hair or jewelry (e.g., earrings), or conductive filaments can be permanently implanted in the ear cartilage to act as antennae using a process similar to ear piercing. The earpiece 10 communicates with a remote processor unit (RPU) 16 via a primary two-way RF link 17, although many other wireless link media (e.g., ultrasonic or infrared) may be used instead. The RPU is typically worn under clothing (e.g., carried in a pocket or purse), but may also be worn in plain sight (e.g., on a belt) if desired. The RPU 16 may be connected (via wired or wireless means 18) to optional secondary wireless link circuitry 19 that allows wireless communication between the RPU and other sources of information (e.g., the general subscriber telephone network) via a secondary wireless link. Note that the optional secondary wireless link circuitry 19 may or may not be contained within the RPU case 70 of FIG. 7.</p> <p>Anderson '783 at 5:3–21</p> <p>During normal operation, speech signals from a nearby talker (and other signals in the ambient audio environment) are picked up by the earpiece wireless microphone 20 and transmitted to the RPU FM receiver 24. The RPU FM receiver 24 output level may be adjusted using the RPU FM receiver 24 volume control. The resulting electrical waveform representing signals from the nearby talker travels through the DPDT switch 25, which is set to the lower position as shown in FIG. 2, to the signal enhancer 26. The signal enhancer 26 may be, for example, a voice changer device that varies the pitch of a received speech signal according to the setting of pushbutton controls located on the RPU signal enhancer 26. The speech signal's pitch can then be raised or lowered as desired to help compensate for a user's hearing loss relative to a particular talker's voice characteristics. The modified speech signal travels from the RPU signal enhancer 26 to the RPU FM transmitter 27, and finally to the earpiece headset 21 receiver where the signal is converted to acoustic waves heard by the user.</p> <p>Anderson '783 at 6:8–25</p> <p>The user's voice (and other signals in the ambient audio environment) is picked up by the earpiece wireless microphone 20, transmitted to the RPU 23 FM receiver 24, sent through the DPDT switch 25 to the microphone input of the handset 28 (the handset's microphone transducer may be removed) and transmitted on the secondary wireless link 29 to, e.g., the general subscriber telephone network. Signals from the secondary wireless link 29 are received by the handset 28, sent from the speaker output of the handset 28 (the handset's speaker may be removed) via the DPDT switch 25 to the signal enhancer 26, FM transmitter 27 and finally to the earpiece headset 21 receiver where the signals are converted to acoustic waves heard by the user. Note that the RPU 23 may also include straight-forward connections to</p>

Claims	Prior Art
	<p>a "personal digital assistant" computer, or voice-operated devices connected to the FM receiver 24 output and not shown in FIG. 2.</p> <p>Anderson '783 at 21:50–62</p> <p>In the simplest case, the ASIC 845 reads data from the A/D 880, performs a multiplication of the data by a constant (read from the default parameter storage 870) to provide amplification, compares the result with maximum allowable values (read from the default parameter storage 870) and writes the result to the D/A 850. The ASIC 845 may perform other processing (e.g., adaptive volume control) as well. Although none of the remote processor and communication features are available in stand-alone mode, which uses pre-selected hearing aid parameter settings from the default parameter storage 870 in the earpiece, the earpiece provides a useful amplified audio path from the microphone to the speaker.</p> <p>Anderson '783 at 22:61–23:35</p> <p>Details of normal system operation: During normal operation of the hearing aid with wireless remote processor, sounds from the ambient environment are picked up by a microphone 12 in the earpiece 10 and transmitted with other information (e.g., data representing the earpiece battery 860 voltage as obtained by an A/D 865 converter and transmitted under control of the earpiece ASIC 845 in the auxiliary bit of many successive replies via the reply data register 825) over a primary two-way wireless link 17 to the RPU 16, where the audio signals are enhanced according to the user's needs by the RPU DSP 948. Signal enhancement is achieved via any of a number of techniques well known in the art that use the capabilities of the general-purpose RPU DSP 948. The enhanced audio signals may be combined with other information (e.g., a synthesized voice, generated within the DSP 948, warning that the aforementioned earpiece battery voltage is low) and transmitted from the RPU 16 over the primary wireless link 17 to the earpiece 10, where they are converted by a speaker 15 to sounds that can only be heard by the user 11. Art auxiliary speaker can be supplied in the RPU 16 as a peripheral 950 device that allows people other than the user to hear the same sounds as those produced by the earpiece speaker 15, if desired (as, e.g., when a person other than the user wishes to use the system's optional cellular telephone capabilities). A different, optional secondary two-way wireless link (shown as part of the RPU in FIG. 9 as secondary wireless link circuitry 944, or connected to the RPU by a wired or other, e.g. infrared, wireless link 18 and shown as optional secondary wireless link circuitry 19 in FIG. 1) can be used for communication between the RPU and a cellular telephone system or other sources of information. An RPU keyboard suitable for covert data entry 72, or voice recognition capabilities in the RPU (implemented, e.g., using any of a number of techniques well known in the art that use the general-purpose RPU DSP 948), can be used to control hearing aid parameters (such as amplification level set by the RPU DSP 948) and telephone dialing</p>

Claims	Prior Art
	<p>functions (e.g., automatic dialing of a predefined telephone number using information stored in the RPU DSP 948) in a manner that can be made imperceptible to a casual observer if desired.</p> <p>Anderson '783 at 24:47–66</p> <p>Morse code keys 956 or voice recognition capabilities implemented using the RPU DSP 948) in a manner that can be made imperceptible to a casual observer if desired. In this respect, the system acts as a "virtual handset" communication aid that replaces a conventional telephone handset in many applications (although users may additionally communicate with a conventional handset or speakerphone in the usual fashion). To accomplish the virtual handset function, the RPU 16 is equipped with optional secondary two-way wireless link circuitry 19 to communicate with a telephone or other communication network. The secondary wireless link is separate from, and does not interfere with, the primary wireless link 17 between the earpiece 10 and RPU 16. In applications where relatively high-power transmissions take place on the secondary link and the secondary link transmitter is integrated into the RPU, the RPU can be placed on a user's wrist or ankle away from vital organs to reduce any potential health risk relative to conventional wireless handsets.</p> <p>Anderson '783 at 25:15–27</p> <p>To obtain solicited information, user requests are processed in a variety of means well known in the art, including inputs via one or more pushbuttons on the RPU and voice recognition computer programs implemented in the RPU electronics. Voice recognition can be used almost exclusively, if desired, thereby eliminating nearly all RPU pushbutton controls. Voice commands recognized by the RPU can cause the RPU to initiate actions such as changing hearing aid parameter settings (e.g., volume), or activate appliances and environmental controls (e.g., aids for the physically handicapped) in a distributed wireless network via the secondary wireless link.</p> <p>A POSITA would have been motivated to combine the teachings of the DHA system with Anderson, to meet the claim limitations with a reasonable expectation of success for the reasons described above with respect to claim 1b.</p> <p>To the extent that the DHA system, alone or in view of Anderson, does not disclose any elements of this claim limitation, it would have been obvious to a person having ordinary skill in the art to combine, with the teachings of DHA systems, those elements, as known to them or as disclosed by other prior art identified in the cover pleading.</p>

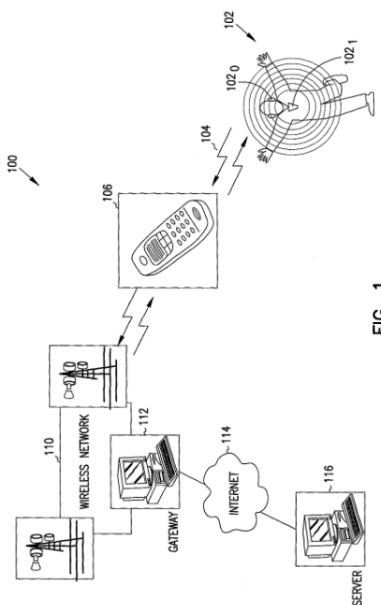
Claims	Prior Art																		
28. A personal audio system as recited in claim 27, wherein the wireless network is a short-range wireless network.	<p>DHA systems, alone or in view of Eaton, expressly or inherently discloses “[a] personal audio system as recited in claim 27, wherein the wireless network is a short-range wireless network.”</p> <p>Phonak SmartLink SX Product Information at p. 2</p> <table><tr><th>Standards</th><th>Europe</th><th>USA</th></tr><tr><td>Wireless microphones</td><td>EN 300.422</td><td>FCC Part 15 Subpart C</td></tr><tr><td>Eletromagn. compatibility</td><td>EN 301.489</td><td>FCC Part 15 Subpart C</td></tr><tr><td>Short-range wireless transmission systems</td><td>EN 300.328</td><td>FCC Part 15 Subpart C</td></tr><tr><td>Safety regulations</td><td>EN 60950</td><td>FCC Part 15 Subpart C</td></tr><tr><td>Bluetooth</td><td colspan="2">*Bluetooth Special Interest Group* conformity test</td></tr></table> <p>Eaton ’458 at Abstract</p> <p>Systems, structures, and methods are provided to fit, program, or upgrade a hearing aid system to a patient. One embodiment includes the use of a mobile device to interact with the hearing aid system through a short-range network. The mobile device is also adapted to communicate with a server through a long-range wireless network. The server may reside on the Internet.</p> <p>Eaton ’458 at 6:11–25</p> <p>The hearing aid system 102 can communicate bi-directionally via a short-range network 104. In one embodiment, the short-range network is a pico-cellular network. In another embodiment, the short-range network includes a network occupying an un-licensed frequency-band. In one embodiment, the short-range network 104 includes wired networking. In another embodiment, the short-range network 104 includes a wireless short-range network such as a radio transmission network or an optical transmission network. One implementation of the radio transmission network includes Bluetooth technology. Bluetooth technology provides a short-range, low-cost radio communication link. Bluetooth may be used to replace wired cables that connect peripherals to this sample of equipment: cellular phones, audio headsets, computer laptops, personal digital assistants, digital cameras, etc. Another implementation of the radio transmission network includes HomeRF, DECT, PHS, or Wireless LAN (WLAN), or other equivalent proprietary wireless communications protocols that do not depart from the present invention.</p>	Standards	Europe	USA	Wireless microphones	EN 300.422	FCC Part 15 Subpart C	Eletromagn. compatibility	EN 301.489	FCC Part 15 Subpart C	Short-range wireless transmission systems	EN 300.328	FCC Part 15 Subpart C	Safety regulations	EN 60950	FCC Part 15 Subpart C	Bluetooth	*Bluetooth Special Interest Group* conformity test	
Standards	Europe	USA																	
Wireless microphones	EN 300.422	FCC Part 15 Subpart C																	
Eletromagn. compatibility	EN 301.489	FCC Part 15 Subpart C																	
Short-range wireless transmission systems	EN 300.328	FCC Part 15 Subpart C																	
Safety regulations	EN 60950	FCC Part 15 Subpart C																	
Bluetooth	*Bluetooth Special Interest Group* conformity test																		

Claims	Prior Art
	<p>Eaton '458 at FIG. 1</p>  <p>A POSITA would have been motivated to combine the teachings of the DHA system with Anderson and Eaton to meet the claim limitations with a reasonable expectation of success for the reasons described above with respect to claim 1b.</p> <p>To the extent that the DHA system, alone or in view of Anderson and Eaton, does not disclose any elements of this claim limitation, it would have been obvious to a person having ordinary skill in the art to combine, with the teachings of DHA systems, those elements, as known to them or as disclosed by other prior art identified in the cover pleading.</p>
<p>29. A personal audio system as recited in claim 27, wherein the short-range wireless network is a Bluetooth network.</p>	<p>DHA systems, alone or in view of Eaton, expressly or inherently discloses “[a] personal audio system as recited in claim 27, wherein the short-range wireless network is a Bluetooth network.”</p>

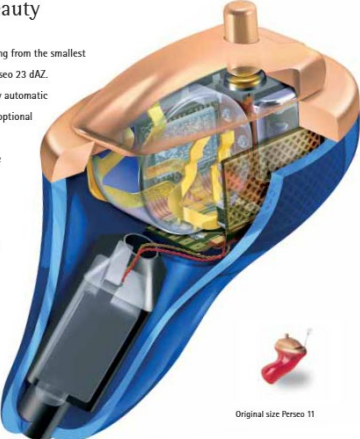
Claims	Prior Art
	<p data-bbox="583 175 1270 211">Phonak - SmartLink SX - FM Systems Manual at p. 3</p> <p data-bbox="638 227 1094 279">SmartLink SX is the first communication device to incorporate</p> <ul data-bbox="638 315 1031 522" style="list-style-type: none"> • Advanced FM features • Hearing Instrument remote control • A Bluetooth link for cell phone use <p data-bbox="583 548 1270 584">Phonak - SmartLink SX - FM Systems Manual at p. 4</p>  <p>The diagram illustrates the SmartLink SX device, a silver, rectangular communication device. It features three circular buttons on its front face, each with a different colored indicator (white, green, and red). Three blue arrows point towards the device from the left, labeled 'Speech', 'Bluetooth link', and 'Audio input'. The 'Speech' arrow originates from a circular inset showing a person's mouth. The 'Bluetooth link' arrow originates from a circular inset showing a mobile phone. The 'Audio input' arrow originates from a circular inset showing a hearing aid. Three horizontal blue lines extend from the right side of the device.</p>

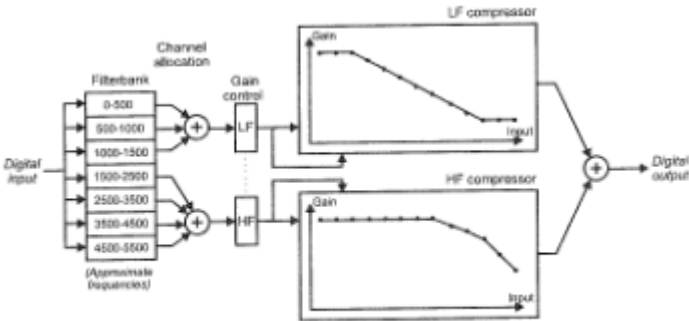
Claims	Prior Art
	<p>Phonak - SmartLink SX - FM Systems Manual at p. 10</p> <div data-bbox="638 224 1602 540"> <p>Smart links with Bluetooth</p> <p>With the inclusion of Bluetooth technology in SmartLink, Phonak now enables hearing instrument wearers to enjoy unrestricted use of all Bluetooth compatible cellular phones. One of the major complaints of hearing impaired people is now history – talking on the telephone is even easier for them than for people with normal hearing. Simply synchronizing the Bluetooth telephone and SmartLink creates the communication path. As long as the cell phone is within a 10-meter range of SmartLink, all incoming calls are automatically directed via SmartLink to the hearing instrument. The cell phone is only needed for dialing outgoing calls. Whenever the cell phone is switched on, the Bluetooth function remains active for incoming calls, regardless of whether SmartLink is switched on or off. This guarantees that telephone calls are not missed, even when the FM function is not required. When using an external audio input to listen to music, incoming telephone calls have priority when the cell phone is switched on. The audio signal is automatically decreased so that the user hears the ringing tone clearly.</p> </div> <p>Phonak - SmartLink SX - FM Systems Manual at p. 10</p>  <p>0</p>

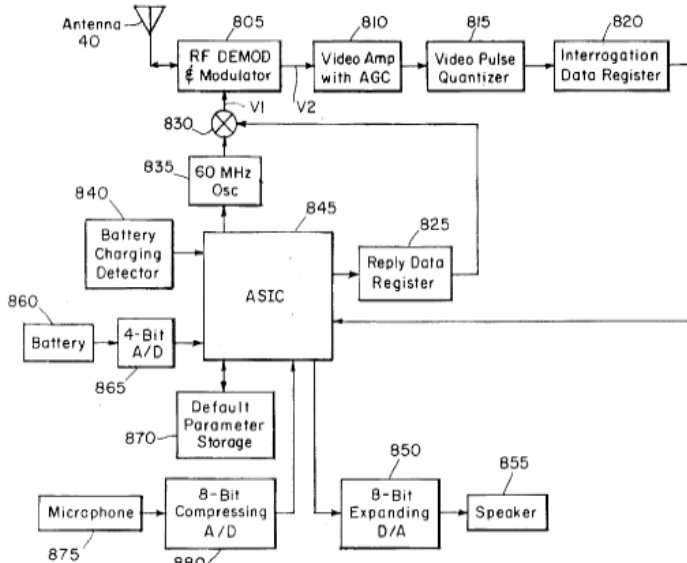
Claims	Prior Art																		
	<p>Phonak - SmartLink SX - FM Systems Manual at p. 14</p> <div></div> <p>Phonak SmartLink SX Product Information at p. 2</p> <table><tr><th>Standards</th><th>Europe</th><th>USA</th></tr><tr><td>Wireless microphones</td><td>EN 300.422</td><td>FCC Part 15 Subpart C</td></tr><tr><td>Eletromagn. compatibility</td><td>EN 301.489</td><td>FCC Part 15 Subpart C</td></tr><tr><td>Short-range wireless transmission systems</td><td>EN 300.328</td><td>FCC Part 15 Subpart C</td></tr><tr><td>Safety regulations</td><td>EN 60950</td><td>FCC Part 15 Subpart C</td></tr><tr><td>Bluetooth</td><td colspan="2">"Bluetooth Special Interest Group" conformity test</td></tr></table> <p>Eaton '458 at 6:11–25</p> <p>The hearing aid system 102 can communicate bi-directionally via a short-range network 104. In one embodiment, the short-range network is a pico-cellular network. In another embodiment, the short-range network includes a network occupying an un-licensed frequency-band. In one embodiment, the short-range network 104 includes wired networking. In another embodiment, the short-range network 104 includes a wireless short-range network such as a radio transmission network or an optical transmission network. One implementation of the radio transmission network includes Bluetooth technology. Bluetooth</p>	Standards	Europe	USA	Wireless microphones	EN 300.422	FCC Part 15 Subpart C	Eletromagn. compatibility	EN 301.489	FCC Part 15 Subpart C	Short-range wireless transmission systems	EN 300.328	FCC Part 15 Subpart C	Safety regulations	EN 60950	FCC Part 15 Subpart C	Bluetooth	"Bluetooth Special Interest Group" conformity test	
Standards	Europe	USA																	
Wireless microphones	EN 300.422	FCC Part 15 Subpart C																	
Eletromagn. compatibility	EN 301.489	FCC Part 15 Subpart C																	
Short-range wireless transmission systems	EN 300.328	FCC Part 15 Subpart C																	
Safety regulations	EN 60950	FCC Part 15 Subpart C																	
Bluetooth	"Bluetooth Special Interest Group" conformity test																		

Claims	Prior Art
	<p>technology provides a short-range, low-cost radio communication link. Bluetooth may be used to replace wired cables that connect peripherals to this sample of equipment: cellular phones, audio headsets, computer laptops, personal digital assistants, digital cameras, etc. Another implementation of the radio transmission network includes HomeRF, DECT, PHS, or Wireless LAN (WLAN), or other equivalent proprietary wireless communications protocols that do not depart from the present invention.</p> <p>Eaton '458 at FIG. 1</p>  <p>A POSITA would have been motivated to combine the teachings of the DHA system with Anderson and Eaton to meet the claim limitations with a reasonable expectation of success for the reasons described above with respect to claim 1b.</p> <p>To the extent that the DHA system, alone or in view of Anderson and Eaton, does not disclose any elements of this claim limitation, it would have been obvious to a person having ordinary skill in the art to combine, with the teachings of DHA systems, those elements, as known to them or as disclosed by other prior art identified in the cover pleading.</p>

Claims	Prior Art
<p>30. A personal audio system as recited in claim 24, wherein the head-worn wearable electronic device includes an accelerometer.</p>	<p>DHA systems, alone or in view of Abreu or Irvin, expressly or inherently discloses “[a] personal audio system as recited in claim 24, wherein the head-worn wearable electronic device includes an accelerometer.”</p> <p>Abreu ’960 at 13:46–60</p> <p>The eyeglasses or other support structures can also have elements which produce and radiate recognizable signals and this procedure could be used to locate and track individuals, particularly in military operations. A permanent magnet can also be mounted in the eyeglasses and used for tracking as described above. A fixed frequency transmitter can be mounted in the eyeglasses and used as a tracking device which utilizes a satellite tracking system by noting the frequency received from the fixed frequency transmitter to a passing satellite, or via Global Positioning Systems. Motion and deceleration can be detected by mounting an accelerometer in the eyeglasses. The use of eyeglasses as tracking means can be useful for locating a kidnapped individual or for rescue operations in the military, since eyeglasses are normally unsuspecting articles.</p> <p>Irvin ’332 at 5:1–23</p> <p>Because the wireless headset 100 of the present invention automatically enables and disables its communication functions based on use, the sensor 220 may be based on a variety of technologies. In general, its purpose is to provide a signal responsive to a natural stimulus associated with use of the headset by a user. Movement or vibration of the wireless headset 100, or placement of the wireless headset 100 upon the user's body indicates usage. Movement may be detected based on vibration or by sensing acceleration, such as with an accelerometer. Mechanical or liquid-based switches represent alternate mechanisms for sensing movement. In such switches, movement causes an electrical connection to be established between contact pairs. Further, the particular sensor 220 employed may be tuned or configured such that only stimulus consistent with actual use actuates the sensor, thus avoiding inadvertent or unintended enabling of the wireless headset 100. In such configurations, the sensor 220 preferably employs output conditioning that prevents assertion of the sensor output signal for stimulus below a defined threshold. Alternatively, the power control circuit 210 may include such conditioning as part of its interface to the sensor 220. FIGS. 4A–4D illustrate possible sensor embodiments.</p> <p>A POSITA would have been motivated to combine the teachings of the DHA system with Abreu or Irvin to meet the claim limitations with a reasonable expectation of success for the reasons described above with respect to claim [1a] and [23a].</p> <p>To the extent that the DHA system, alone or in view of Abreu or Irvin, does not disclose any elements of this claim limitation, it would have been obvious to a person having ordinary skill in the art to combine, with the</p>

Claims	Prior Art
	<p>teachings of DHA systems, those elements, as known to them or as disclosed by other prior art identified in the cover pleading.</p>
<p>31. A personal audio system as recited in claim 24, wherein the head-worn wearable electronic device includes a plurality of circuit boards, each circuit board including at least one integrated circuit mounted thereto.</p>	<p>DHA systems, alone or in view of Anderson and Jannard, expressly or inherently discloses “[a] personal audio system as recited in claim 24, wherein the head-worn wearable electronic device includes a plurality of circuit boards, each circuit board including at least one integrated circuit mounted thereto.”</p> <p>Phonak – Perseo Brochure at Page 22</p> <div data-bbox="646 448 1276 1045"> <h2 style="text-align: center;">The beauty of high technology</h2> <p style="text-align: center;">Perseo[®] ITEs – hidden beauty</p> <div style="display: flex; align-items: center;"> <div style="flex: 1;"> <p>The Perseo family includes a full line of ITEs ranging from the smallest Perseo 11 CIC to the elegant dual-microphone Perseo 23 dAZ.</p> <p>All CIC and ITE instruments, like the BTEs, are fully automatic and have PersonalLogic. A TacTronic switch is optional on all ITE models and easily added after fitting by exchanging battery covers. This is a very simple and practical means of providing extra control and additional programs.</p> <p>In many countries, ITEs are available with Phonak eShell. These electronically created shells are hypoallergenic, exceptionally durable and more easily retained in the wearer's ear.</p> </div> <div style="flex: 1; text-align: center;">  </div> </div> </div> <p>Oticon – DigiFocus Compact User Instructional Brochure at p. 2</p> <p style="text-align: center;">The programming is accomplished via computer software (OtiSet) using the NOAH standard platform. Additionally, DigiSound can be programmed from Oticon's dedicated programming box, “Easy fit”.</p>

Claims	Prior Art
	<p>Oticon – DigiFocus Description at Page 1</p> <p>Inside each DigiFocus, the digital chips have the same processing power as a fast laptop computer. This power is used to analyse the incoming sound signal, and speech is extracted through a unique system called ASA - Adaptive Speech Alignment. ASA always makes speech audible, amplifying it as much as necessary without making it uncomfortably loud or distorted. Sounds are reproduced very naturally - making speech far easier to understand than with any other hearing instrument on the market.</p> <p>Oticon – DigiFocus White Paper (1996) at Page 1</p> <p>The OtiSet software is the single most important tool for fitting DigiFocus. Underneath its simple and easy to use interface, OtiSet sums up all of the advanced knowledge on which DigiFocus is built. By analyzing all of the available personal data, OtiSet automatically recommends the most appropriate settings. In short, OtiSet enables you to draw on all the audiological knowledge available from Oticon. □</p> <p>Oticon – DigiFocus White Paper (1996) at Fig. 1</p>  <p><i>Fig.1. DigiFocus - block diagram of the digital audio processor, DAP.</i></p> <p>Oticon – LA Time Article re THE GOODS: The Sound Wave of the Future (Nov. 10 1995) at 1</p> <p>DigiFocus operates on two tiny new computer chips that have the processing power of a normal desktop computer. Traditional hearing instruments have only up to 10 adjustable parameters, but Oticon representatives say that their processor will give a maximum of 100 software-controlled adjustment parameters. And the device can be programmed by an audiologist while the hearing-impaired person is wearing it.</p>

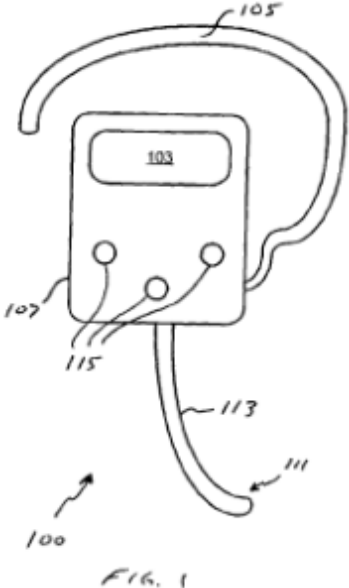
Claims	Prior Art
	<p data-bbox="583 175 898 207">Anderson '783 at FIG. 8</p>  <p data-bbox="940 795 1018 824">FIG. 8</p> <p data-bbox="583 844 919 876">Anderson '783 at 3:40–41</p> <p data-bbox="625 885 1963 950">FIG. 8 is a block diagram showing details of an earpiece for a full-featured preferred embodiment of the invention.</p> <p data-bbox="583 1006 940 1039">Anderson '783 at 12:38–46</p> <p data-bbox="625 1047 1942 1226">Note that all earpiece functions (e.g., power control for the 60 MHz oscillator 835) are controlled by a custom application-specific integrated circuit (ASIC) 845. The ASIC 845 may be a low-power low-voltage device similar, e.g., to the ECI Semiconductor P576 custom programmable linear/digital gate array. Many of the functions that are shown as separate blocks for clarity in FIG. 8, e.g., reply data register 825, may be implemented within the ASIC 845.</p>

Claims	Prior Art
	<p data-bbox="583 175 884 207">Jannard '526 at FIG. 3J</p> <p data-bbox="919 784 1020 816"><i>FIG. 3J</i></p> <p data-bbox="583 846 894 878">Jannard '526 at 3:56–58</p> <p data-bbox="632 883 1990 948">FIG. 3J is a front, top, and left side perspective and exploded view of the wearable audio device illustrated in FIG. 3D.</p> <p data-bbox="583 1008 936 1040">Jannard '526 at 12:62–13:6</p> <p data-bbox="632 1045 1990 1295">A cover member 262 cooperates with the housing 250 to define the closed internal cavity 252. In the illustrated embodiment, the internal cavity 252 includes at least one compartment configured to receive an electronic circuit board 264 which includes at least one switch for each of the buttons 73a, 73b. In an exemplary but non-limiting embodiment, the board 264 can include two switches, one for each of the buttons 73a, 73b, which are configured to control a volume output from the speakers 14A', 16A'. The cover 262 can be attached to the ear stems 54A' with any type of fastener, such as, for example, but without limitation, screws, rivets, bolts, adhesive, and the like.</p> <p data-bbox="583 1356 911 1388">Jannard '526 at 13:50–55</p> <p data-bbox="632 1393 1990 1458">The internal cavity 282, in the illustrated embodiment, is configured to receive a printed circuit board 284. In the illustrated embodiment, the printed circuit board 284 includes one switch for each of the buttons</p>

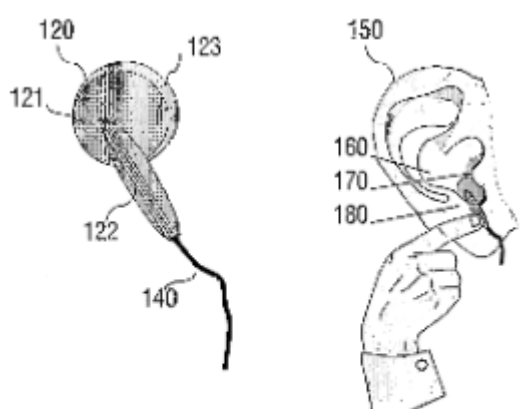
Claims	Prior Art
	<p>73c, 73d, and 73e. Additionally, the printed circuit board 284 includes an audio file storage and playback device 286.</p> <p>Jannard '526 at 14:6–14 Advantageously, the printed circuit board 284 also includes or is in electrical communication with a data transfer port 388. In the illustrated embodiment, the housing 280 includes an aperture (not shown) disposed in a position similar to the position of the aperture 272 on the housing 250. In the housing 280, however, the aperture is aligned with the data transfer port 288. Thus, when the printed circuit board 284 is received in the internal cavity 282, the data transfer port 288 is aligned with the aperture.</p> <p>Jannard '526 at 15:60–16:2 Thus, with reference to FIGS. 3J and 3M, the ribbon connector 318 can extend from the housing 280, into the recesses 294, 296, through an aperture (not shown) in the hinge pin 298 to the upper opening within the hinge pin 298, then through the channel 316 (FIG. 3E), to an upper opening of the hinge pin 270, out through an aperture (not shown) through a side of a hinge pin 270, through the recesses 266, 268 of the housing 250, and then to the speaker 14A', printed circuit board 264, and the power storage device 274.</p> <p>Jannard '526 at 25:12–21 With reference to FIGS. 8, 9A, and 9B, in another embodiment, the audio device 10C is advantageously adapted to support any of a variety of portable electronic circuitry or devices which have previously been difficult to incorporate into conventional headsets due to bulk, weight or other considerations. For example, but without limitation, the electronics are digital or other storage devices and retrieval circuitry such as for retrieving music or other information from MP3 format memory or other memory devices.</p> <p>A POSITA would have been motivated to combine the teachings of the DHA system with Anderson and Jannard to meet the claim limitations with a reasonable expectation of success for the reasons described above with respect to claim [1a] and [26a].</p> <p>To the extent that the DHA system, alone or in view of Anderson and Abreu, does not disclose any elements of this claim limitation, it would have been obvious to a person having ordinary skill in the art to combine, with the teachings of DHA systems, those elements, as known to them or as disclosed by other prior art identified in the cover pleading.</p>

Claims	Prior Art
<p>32. A personal audio system as recited in claim 31, wherein at least one of the circuit boards is a flexible printed circuit board.</p>	<p>DHA systems, alone or in view of Jannard, expressly or inherently discloses “[a] personal audio system as recited in claim 31, wherein at least one of the circuit boards is a flexible printed circuit board.”</p> <p>Jannard ’526 at FIG. 3J</p> <p>FIG. 3J</p> <p>Jannard ’526 at 3:56–58</p> <p>FIG. 3J is a front, top, and left side perspective and exploded view of the wearable audio device illustrated in FIG. 3D.</p> <p>Jannard ’526 at 12:62–13:6</p> <p>A cover member 262 cooperates with the housing 250 to define the closed internal cavity 252. In the illustrated embodiment, the internal cavity 252 includes at least one compartment configured to receive an electronic circuit board 264 which includes at least one switch for each of the buttons 73a, 73b. In an exemplary but non-limiting embodiment, the board 264 can include two switches, one for each of the buttons 73a, 73b, which are configured to control a volume output from the speakers 14A’, 16A’. The cover 262 can be attached to the ear stems 54A’ with any type of fastener, such as, for example, but without limitation, screws, rivets, bolts, adhesive, and the like.</p>

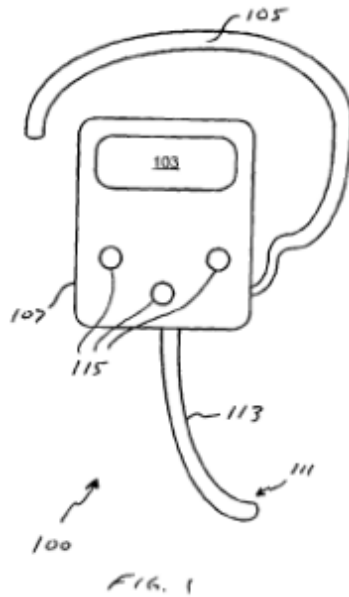
Claims	Prior Art
	<p>Jannard '526 at 14:6–14</p> <p>Advantageously, the printed circuit board 284 also includes or is in electrical communication with a data transfer port 388. In the illustrated embodiment, the housing 280 includes an aperture (not shown) disposed in a position similar to the position of the aperture 272 on the housing 250. In the housing 280, however, the aperture is aligned with the data transfer port 288. Thus, when the printed circuit board 284 is received in the internal cavity 282, the data transfer port 288 is aligned with the aperture.</p> <p>Jannard '526 at 26:21–27</p> <p>With reference to FIGS. 8, 9A, and 9B, in another embodiment, the transceiver 114 is adapted to employ a wide variety of technologies, including wireless communication such as RF, IR, ultrasonic, laser or optical, as well as wired and other communications technologies. In one embodiment, a body-LAN radio is employed. Other embodiments can employ a flexible-circuit design.</p> <p>A POSITA would have been motivated to combine the teachings of the DHA system with Jannard to meet the claim limitations with a reasonable expectation of success for the reasons described above with respect to claim [1a] and [26a].</p> <p>To the extent that the DHA system, alone or in view of Jannard, does not disclose any elements of this claim limitation, it would have been obvious to a person having ordinary skill in the art to combine, with the teachings of DHA systems, those elements, as known to them or as disclosed by other prior art identified in the cover pleading.</p>
<p>33. A personal audio system as recited in claim 24, wherein the head-worn wearable electronic device includes a touch sensitive surface.</p>	<p>DHA systems, alone or in view of Bodley or Hollemans or Seshadri, expressly or inherently discloses “[a] personal audio system as recited in claim 24, wherein the head-worn wearable electronic device includes a touch sensitive surface.”</p>

Claims	Prior Art
	<p data-bbox="583 175 865 207">Bodley '433 at FIG. 1</p>  <p data-bbox="583 846 856 878">Bodley '433 at 2:7–8</p> <p data-bbox="632 883 1915 915">FIG. 1 is a frontal view of a generic headset that includes a display in accordance with the invention;</p> <p data-bbox="583 971 888 1003">Bodley '433 at 2:38–50</p> <p data-bbox="632 1008 1984 1295">Although headset 100 may use any of a variety of means of attaching to a user (e.g., headband, earpiece, etc.) and thus allowing hands-free operation of the headset, preferably headset 100 includes an earpiece 105 (e.g., an ear hook) that is coupled to the user's ear, or a portion thereof. More preferably, earpiece 105 fits around a portion of the user's ear. Earpiece 105 is coupled to body 107 of the headset, body 107 housing the various electronic components necessary to provide functionality to headset 100. Attached to the back surface of body 107 is an output transducer 109 (e.g., a speaker). Also coupled to body 107 is an input transducer 111 (e.g., a microphone), which, in this embodiment, is located at the end of a boom member 113.</p> <p data-bbox="583 1351 913 1383">Bodley '433 at 2:65–3:22</p> <p data-bbox="632 1388 1976 1458">Display 103 can be used to provide the user with various types of information depending upon the design goals for the specific headset. The design goals are primarily based on the expected headset use, desired</p>

Claims	Prior Art
	<p>headset functionality, display size, display resolution, allowable headset power drain, and the processing capabilities of the headset on-board processor (also referred to as a controller or an interface controller). Examples of information that can be displayed on display 103 are headset battery level, headset signal level, call status (e.g., dialing, in use, etc.), caller identification, current time/date, and timer information such as elapsed time associated with a particular call. Display 103 can also be configured to flash or otherwise visually indicate an incoming call. Assuming that the information is provided to the headset processor, display 103 can also display the battery level of the base unit (e.g., cell phone, dongle, etc.) and the signal strength of the base unit (e.g., cell phone). Additionally, display 103 can be used during headset configuration (e.g., to adjust or set the ring mode, ring loudness, headset volume, display mode, voice dialing, time, date, etc.) or to provide additional headset functionality (e.g., phone lists, text messages, dialing using virtual numeric keypad or phone list dialing). Additionally, display 103 can be touch sensitive, thus providing a touch screen for data input (e.g., phone numbers into a phone list, dialing using a virtual numeric keypad).</p> <p>Hollemans '490 at 2:10–18</p> <p>The remote controller thus gets a fixed position with respect to the user in or by the human ear. It is easy to find one's ear. This prevents looking for the remote controller. In addition, the user is not tempted to look at the remote controller, because there is nothing relevant to be seen for controlling the device. To avoid having to find tiny buttons, merely touching the remote controller somewhere on its touch-sensitive area operates said remote controller, which detects being touched and subsequently sends the control signal to the device. Also the problem of finding the proper orientation of the remote controller is solved by wearing the remote controller in or by the ear, because its orientation becomes fixed with respect to the user.</p> <p>Hollemans '490 at 2:25–33</p> <p>Advantageously, the controller is arranged to detect a temporal pattern in the touch-sensitive area being touched, and to send the control signal in response to detecting the temporal pattern. In general, the device offers a plurality of functions and capabilities. Rather than having more buttons on the remote control that may be hard to distinguish, the single touch-sensitive area is used for controlling the plurality of functions. The temporal pattern is a particularly appropriate user interface, because it is easy to create temporal patterns in touching the area rhythmically and because the temporal patterns can constitute a natural and consistent interface. An example is a single short tap for pause/play, double short tap for next track, triple short tap for next artist/album and a long tap for adjusting the volume.</p>

Claims	Prior Art
	<p>Hollemans '490 at 3:1–6</p> <p>Advantageously, the outer surface has a further touch-sensitive area, such that the further touch-sensitive area is touched substantially by the ear when the controller is substantially worn in or by a human ear, the controller being arranged to send the control signal only if the further touch-sensitive area is touched. This will prevent the accidental remote control of the device when the remote controller is not worn, thus preventing unwanted activation of the device and unnecessary battery exhaustion.</p> <p>Hollemans '490 at FIG. 2</p>  <p>FIG.2</p> <p>Seshadri '839 at [0041]</p> <p>Wearable earpiece 12, as shown in FIG. 5, includes antenna 46, transmit/receive switch 122, receiver module 40, data recovery module 42, speaker module 44, transmitter module 120, input module 128 and display module 132. Receiver module 40, data recovery module 42 and speaker module 44 operate as discussed with reference to FIG. 4, with the exception of the data recovery module 42, to produce display information that is provided to display module 132. For instance, the received RF signal may include display information such as caller ID, command information, etc. which is separated by data recovery module 42 and provided to display module 132, which may be an LCD display, plasma display, etc.</p>

Claims	Prior Art
	<p>Seshadri '839 at [0042] Input module 128, which may be a keypad, touch screen, voice recognition circuit, or other like user interfaces, receives user commands and produces digital command messages 124 there from. Such digital command messages 124 includes, but are not limited to, packet size, synchronization information, frequency hopping initiation information, timeslot allocation information, link establishment information, piconet address information, fast-forward, play, pause, volume adjust, record, stop and rewind.</p> <p>A POSITA would have been motivated to combine the teachings of the DHA system and Anderson and Bodley or Hollemans or Seshadri to meet the claim limitations with a reasonable expectation of success for the reasons described above with respect to claim [1a]. The DHA system discloses that a user can control various functions and settings of an earpiece for a hearing aid system. A POSITA would recognize that Anderson discloses similar features, and also that the earpiece can occasionally lose connection to the RPU, in which case it can operate in a stand-alone mode. When the earpiece is operating in standalone mode, the earpiece settings (e.g., hearing aid parameters) can only be adjusted when the earpiece is being charged. Thus, Anderson '783's hearing aid system faces the problem of how to adjust operation of the earpiece without the RPU while the earpiece is in use. A POSITA would recognize Bodley '433 or Seshadri to provide a solution to this problem: providing the earpiece with a touch sensitive surface that can be used for controlling functions and settings of the earpiece. Moreover, the DHA system, Anderson '783 and Bodley '433 and Seshadri share various structural and functional similarities that render the inventions readily combinable. For example, each discloses ear-worn devices for providing audio output to a user. Each includes similar components, such as microphones, speakers, wireless transceivers, and control circuitry. Accordingly, combining the teachings of the DHA system with Anderson '783 and Bodley '433 or Seshadri would involve combining known components to yield predictable results.</p> <p>To the extent that the DHA system, alone or in view of Anderson and Bodley or Seshadri, does not disclose any elements of this claim limitation, it would have been obvious to a person having ordinary skill in the art to combine, with the teachings of DHA systems, those elements, as known to them or as disclosed by other prior art identified in the cover pleading.</p>
34. A personal audio system as recited in claim 33, wherein the touch sensitive surface is configured to detect a user input, the detected user input operating to	DHA systems, alone or in view of Bodley or Hollemans or Seshadri, expressly or inherently discloses “[a] personal audio system as recited in claim 33, wherein the touch sensitive surface is configured to detect a user input, the detected user input operating to at least alter operation of the head-worn wearable electronic device.”

Claims	Prior Art
<p>at least alter operation of the head-worn wearable electronic device.</p>	<p>Bodley '433 at FIG. 1</p>  <p>Bodley '433 at 2:7–8</p> <p>FIG. 1 is a frontal view of a generic headset that includes a display in accordance with the invention;</p> <p>Bodley '433 at 2:38–50</p> <p>Although headset 100 may use any of a variety of means of attaching to a user (e.g., headband, earpiece, etc.) and thus allowing hands-free operation of the headset, preferably headset 100 includes an earpiece 105 (e.g., an ear hook) that is coupled to the user's ear, or a portion thereof. More preferably, earpiece 105 fits around a portion of the user's ear. Earpiece 105 is coupled to body 107 of the headset, body 107 housing the various electronic components necessary to provide functionality to headset 100. Attached to the back surface of body 107 is an output transducer 109 (e.g., a speaker). Also coupled to body 107 is an input transducer 111 (e.g., a microphone), which, in this embodiment, is located at the end of a boom member 113.</p> <p>Bodley '433 at 2:65–3:22</p> <p>Display 103 can be used to provide the user with various types of information depending upon the design goals for the specific headset. The design goals are primarily based on the expected headset use, desired</p>

Claims	Prior Art
	<p>headset functionality, display size, display resolution, allowable headset power drain, and the processing capabilities of the headset on-board processor (also referred to as a controller or an interface controller). Examples of information that can be displayed on display 103 are headset battery level, headset signal level, call status (e.g., dialing, in use, etc.), caller identification, current time/date, and timer information such as elapsed time associated with a particular call. Display 103 can also be configured to flash or otherwise visually indicate an incoming call. Assuming that the information is provided to the headset processor, display 103 can also display the battery level of the base unit (e.g., cell phone, dongle, etc.) and the signal strength of the base unit (e.g., cell phone). Additionally, display 103 can be used during headset configuration (e.g., to adjust or set the ring mode, ring loudness, headset volume, display mode, voice dialing, time, date, etc.) or to provide additional headset functionality (e.g., phone lists, text messages, dialing using virtual numeric keypad or phone list dialing). Additionally, display 103 can be touch sensitive, thus providing a touch screen for data input (e.g., phone numbers into a phone list, dialing using a virtual numeric keypad).</p> <p>Seshadri '839 at [0041] Wearable earpiece 12, as shown in FIG. 5, includes antenna 46, transmit/receive switch 122, receiver module 40, data recovery module 42, speaker module 44, transmitter module 120, input module 128 and display module 132. Receiver module 40, data recovery module 42 and speaker module 44 operate as discussed with reference to FIG. 4, with the exception of the data recovery module 42, to produce display information that is provided to display module 132. For instance, the received RF signal may include display information such as caller ID, command information, etc. which is separated by data recovery module 42 and provided to display module 132, which may be an LCD display, plasma display, etc.</p> <p>Seshadri '839 at [0042] Input module 128, which may be a keypad, touch screen, voice recognition circuit, or other like user interfaces, receives user commands and produces digital command messages 124 there from. Such digital command messages 124 includes, but are not limited to, packet size, synchronization information, frequency hopping initiation information, timeslot allocation information, link establishment information, piconet address information, fast-forward, play, pause, volume adjust, record, stop and rewind.</p>

Claims	Prior Art
	<p>A POSITA would have been motivated to combine the teachings of the DHA system with Anderson and Bodley or Seshadri to meet the claim limitations with a reasonable expectation of success for the reasons described above with respect to claim [1a] and [33].</p> <p>To the extent that the DHA system, alone or in view of Anderson and Bodley or Hollemans or Seshadri, does not disclose any elements of this claim limitation, it would have been obvious to a person having ordinary skill in the art to combine, with the teachings of DHA systems, those elements, as known to them or as disclosed by other prior art identified in the cover pleading.</p>
<p>[35a] A personal audio system as recited in claim 34, wherein the head-worn wearable electronic device is configured to support audio output to the user, and</p>	<p>DHA systems, alone or in view of Anderson, expressly or inherently discloses “[a] personal audio system as recited in claim 34, wherein the head-worn wearable electronic device is configured to support audio output to the user.”</p> <p><i>See above at claim [1g].</i></p> <p>A POSITA would have been motivated to combine the teachings of the DHA system with Anderson to meet the claim limitations with a reasonable expectation of success for the reasons described above with respect to claim [1g].</p> <p>To the extent that the DHA system, alone or in view of Anderson and Bodley, does not disclose any elements of this claim limitation, it would have been obvious to a person having ordinary skill in the art to combine, with the teachings of DHA systems, those elements, as known to them or as disclosed by other prior art identified in the cover pleading.</p>
<p>[35b] wherein the at least one operation of the head-worn electronic device being altered concerns control of the audio output being provided to the user.</p>	<p>DHA systems, alone or in view of Bodley or Hollemans ’490, expressly or inherently discloses “wherein the at least one operation of the head-worn electronic device being altered concerns control of the audio output being provided to the user.”</p> <p>Phonak – Fitting Guide for MLxi and MLxi Baha at Page 3</p> <p><small>Most Phonak HIs offer the EasyFM feature. If this is enabled, the HI automatically switches to FM mode a few seconds after a signal (e.g. voice or music) has been broadcast by the FM transmitter. The reaction time of the EasyFM feature ranges up to several seconds depending on the hearing instrument and the signal intensity. The HI automatically switches back to the default listening program after approximately one minute, when no signal has been broadcasted by the transmitter.</small></p> <p>Bodley ’433 at 2:65–3:22</p> <p>Display 103 can be used to provide the user with various types of information depending upon the design goals for the specific headset. The design goals are primarily based on the expected headset use, desired headset functionality, display size, display resolution, allowable headset power drain, and the processing</p>

Claims	Prior Art
	<p>capabilities of the headset on-board processor (also referred to as a controller or an interface controller). Examples of information that can be displayed on display 103 are headset battery level, headset signal level, call status (e.g., dialing, in use, etc.), caller identification, current time/date, and timer information such as elapsed time associated with a particular call. Display 103 can also be configured to flash or otherwise visually indicate an incoming call. Assuming that the information is provided to the headset processor, display 103 can also display the battery level of the base unit (e.g., cell phone, dongle, etc.) and the signal strength of the base unit (e.g., cell phone). Additionally, display 103 can be used during headset configuration (e.g., to adjust or set the ring mode, ring loudness, headset volume, display mode, voice dialing, time, date, etc.) or to provide additional headset functionality (e.g., phone lists, text messages, dialing using virtual numeric keypad or phone list dialing). Additionally, display 103 can be touch sensitive, thus providing a touch screen for data input (e.g., phone numbers into a phone list, dialing using a virtual numeric keypad).</p> <p>Hollemans '490 at 2:10–18</p> <p>The remote controller thus gets a fixed position with respect to the user in or by the human ear. It is easy to find one's ear. This prevents looking for the remote controller. In addition, the user is not tempted to look at the remote controller, because there is nothing relevant to be seen for controlling the device. To avoid having to find tiny buttons, merely touching the remote controller somewhere on its touch-sensitive area operates said remote controller, which detects being touched and subsequently sends the control signal to the device. Also the problem of finding the proper orientation of the remote controller is solved by wearing the remote controller in or by the ear, because its orientation becomes fixed with respect to the user.</p> <p>Hollemans '490 at 2:25–33</p> <p>Advantageously, the controller is arranged to detect a temporal pattern in the touch-sensitive area being touched, and to send the control signal in response to detecting the temporal pattern. In general, the device offers a plurality of functions and capabilities. Rather than having more buttons on the remote control that may be hard to distinguish, the single touch-sensitive area is used for controlling the plurality of functions. The temporal pattern is a particularly appropriate user interface, because it is easy to create temporal patterns in touching the area rhythmically and because the temporal patterns can constitute a natural and consistent interface. An example is a single short tap for pause/play, double short tap for next track, triple short tap for next artist/album and a long tap for adjusting the volume.</p>

Claims	Prior Art
	<p>Hollemans '490 at 3:1–6</p> <p>Advantageously, the outer surface has a further touch-sensitive area, such that the further touch-sensitive area is touched substantially by the ear when the controller is substantially worn in or by a human ear, the controller being arranged to send the control signal only if the further touch-sensitive area is touched. This will prevent the accidental remote control of the device when the remote controller is not worn, thus preventing unwanted activation of the device and unnecessary battery exhaustion.</p> <p>Hollemans '490 at 3:7–10</p> <p>Advantageously, the controller is arranged to send a further control signal to the device if the further touch-sensitive area is touched. This makes it, for example, possible for the device to pause playing music as soon as the earpiece is taken out and to resume playing music upon re-insertion.</p> <p>Hollemans '490 at 3:11–21</p> <p>In an embodiment of the personal audio system according to the invention, the system comprises a second controller for remotely controlling the device by sending a further control signal to the device, the second controller having an outer surface with a further touch-sensitive area, the second controller being arranged to be substantially worn in or by a human ear, and the second controller being further arranged to detect a further temporal pattern in the further touch-sensitive area being touched, and to send the further control signal in response to detecting the further temporal pattern. A second controller fits in with the natural symmetry of the human head. It also considerably enhances the user interface for controlling the device. For example, the symmetry can be exploited in such a way that a long tap at the left remote controller decreases the volume, but a long tap at the right controller increases the volume.</p> <p>Hollemans '490 at FIG. 3</p>  <p style="text-align: center;">FIG.3</p>

Claims	Prior Art
	<p>Holleman '490 at FIG. 4</p>  <p style="text-align: center;">FIG.4</p> <p>A POSITA would have been motivated to combine the teachings of the DHA system with Anderson and Bodley or Hollemans to meet the claim limitations with a reasonable expectation of success for the reasons described above with respect to claim [1a] and [33].</p> <p>To the extent that the DHA system, alone or in view of Anderson and Bodley or Hollemans, does not disclose any elements of this claim limitation, it would have been obvious to a person having ordinary skill in the art to combine, with the teachings of DHA systems, those elements, as known to them or as disclosed by other prior art identified in the cover pleading.</p>
<p>36. A personal audio system as recited in claim 24, wherein the head-worn wearable electronic device includes a force sensor configured to receive a user input based at least on sensing force applied to the force sensor.</p>	<p>DHA systems, alone or in view of Abreu or Irvin, expressly or inherently discloses “[a] personal audio system as recited in claim 24, wherein the head-worn wearable electronic device includes a force sensor configured to receive a user input based at least on sensing force applied to the force sensor.”</p> <p>Abreu '960 at 15:60–65</p> <p>By way of illustration, but not of limitation, pressure sensing means can be incorporated in the end of the temples to detect if the sunglasses are being worn, and an abrupt drop in the pressure signal indicates glasses were removed or misplacement of the sensor can also generate an identifiable signal.</p> <p>Irvin '332 at 6:34–48</p> <p>FIG. 4E illustrates the power control circuit 210 coupled to an alternate type of attitude sensor 220. In this embodiment, the attitude sensor 220 is a mechanical sensor that includes a contacting member 224 responsive to a physical orientation of the sensor 220. In one or more physical orientations, the contacting member 224 does not electrically connect the first and second contacts of attitude sensor 220. In one or more other physical orientations, the contacting member places the first and second contacts of attitude sensor 220 in electrical connection. The power control circuit 210 is configured to be responsive to this</p>

Claims	Prior Art
	<p>electrical connection. Contacting member 224 may be resiliently biased away from one or both contacts to ensure an open electrical connection between contacts when the headset is not oriented in a position indicative of use.</p> <p>A POSITA would have been motivated to combine the teachings of the DHA system with Abreu or Irvin to meet the claim limitations with a reasonable expectation of success for the reasons described above with respect to claim [1a] and [23a].</p> <p>To the extent that the DHA system, alone or in view of Abreu or Irvin, does not disclose any elements of this claim limitation, it would have been obvious to a person having ordinary skill in the art to combine, with the teachings of DHA systems, those elements, as known to them or as disclosed by other prior art identified in the cover pleading.</p>
<p>37. A personal audio system as recited in claim 36, wherein the force sensor is configured to detect the user input, the detected user input operating to alter at least an operation of the head-worn wearable electronic device.</p>	<p>DHA systems, alone or in view of Abreu or Irvin, expressly or inherently discloses “[a] personal audio system as recited in claim 36, wherein the force sensor is configured to detect the user input, the detected user input operating to alter at least an operation of the head-worn wearable electronic device.”</p> <p>Abreu ’960 at 15:60–65</p> <p>By way of illustration, but not of limitation, pressuse sensing means can be incorporated in the end of the temples to detect if the sunglasses are being worn, and an abrupt drop in the pressure signal indicates glasses were removed or misplacement of the sensor can also generate an identifiable signal.</p> <p>Irvin ’332 at 2:20–34</p> <p>The wireless headset serves as a communication accessory to a base unit, such as a cellular telephone. When not in use, the wireless headset reverts to an inactive state and draws essentially no power. A power control circuit imparts automatic enabling capability to the wireless headset. Preferably, the power control circuit employs a sensor responsive to a natural stimulus associated with a person preparing to use the wireless headset. The wireless headset may sense motion, physical orientation, or user proximity. Sensing techniques include but are not limited to micro-electro-mechanical-systems (MEMS), mercury switches, attitude switches, and IR or acoustic proximity sensors. In some embodiments, the wireless headset transitions from the inactive state to an active state in response to a user picking it up or mounting the wireless headset upon their head.</p> <p>Irvin ’332 at 6:34–48</p> <p>FIG. 4E illustrates the power control circuit 210 coupled to an alternate type of attitude sensor 220. In this</p>

Claims	Prior Art
	<p>embodiment, the attitude sensor 220 is a mechanical sensor that includes a contacting member 224 responsive to a physical orientation of the sensor 220. In one or more physical orientations, the contacting member 224 does not electrically connect the first and second contacts of attitude sensor 220. In one or more other physical orientations, the contacting member places the first and second contacts of attitude sensor 220 in electrical connection. The power control circuit 210 is configured to be responsive to this electrical connection. Contacting member 224 may be resiliently biased away from one or both contacts to ensure an open electrical connection between contacts when the headset is not oriented in a position indicative of use.</p> <p>Irvin '332 at 8:11–22</p> <p>FIGS. 6 and 7 illustrate operation of the wireless headset 100. FIG. 6 is a state diagram for one exemplary embodiment of the present invention. When the wireless headset 100 is not in use, it automatically enters an inactive state (610). As explained, this inactive state may constitute disconnecting the communications circuit 230 from the power source 240, or may simply represent placing the communications circuit 230 in a non-functional, low-power condition. Absent user stimulus, such as motion of the wireless headset 100 or placement of the wireless headset 100 upon the user's body, the wireless headset 100 remains in the inactive state (610).</p> <p>A POSITA would have been motivated to combine the teachings of the DHA system with Abreu or Irvin to meet the claim limitations with a reasonable expectation of success for the reasons described above with respect to claim [1a] and [23a].</p> <p>To the extent that the DHA system, alone or in view of Abreu or Irvin, does not disclose any elements of this claim limitation, it would have been obvious to a person having ordinary skill in the art to combine, with the teachings of DHA systems, those elements, as known to them or as disclosed by other prior art identified in the cover pleading.</p>
<p>38[a] A personal audio system as recited in claim 37, wherein the head-worn wearable electronic device is configured to support audio output to the user, and</p>	<p>DHA systems, alone or in view of Anderson, expressly or inherently discloses “[a] personal audio system as recited in claim 37, wherein the head-worn wearable electronic device is configured to support audio output to the user.”</p> <p><i>See</i> above at claim [1g].</p> <p>A POSITA would have been motivated to combine the teachings of the DHA system with Anderson and Abreu to meet the claim limitations with a reasonable expectation of success for the reasons described above with respect to claim [1g] and [23a].</p>

Claims	Prior Art
	<p>To the extent that the DHA system, alone or in view of Anderson, does not disclose any elements of this claim limitation, it would have been obvious to a person having ordinary skill in the art to combine, with the teachings of DHA systems, those elements, as known to them or as disclosed by other prior art identified in the cover pleading.</p>
<p>[38b] wherein the at least one operation of the head-worn electronic device being altered concerns control of the audio output being provided to the user.</p>	<p>DHA systems, alone or in view of Anderson and Abreu or Irvin, expressly or inherently discloses “wherein the at least one operation of the head-worn electronic device being altered concerns control of the audio output being provided to the user.”</p> <p><i>See</i> above at claim [35b], [23b].</p> <p>Abreu ’960 discloses that the eyeglasses can “send a signal if the eyeglasses were removed,” which can be detected based on “an abrupt drop in the pressure signal.” <i>See id.</i> at 15:57-63. Abreu ’960 further discloses that signals produced by the apparatus can be reported “by audio, visual or tactile means.” <i>See id.</i> at 19:25-26. Thus, the operation of Abreu ’960’s eyeglasses that is altered based on a detected user input to a force sensor can be the generation or cessation of an audio signal. A POSITA would have understood this to “concern” control of the audio output being provided to the user in that it dictates whether or not audio is being output.</p> <p>Similarly, Irvin discloses that entering an “inactive state” upon determining that the device has been removed from the user’s body, and that the device is “non-functional” in the low-power state. Irvin ’332 at 8:11–22. A POSITA would understand that that the operation being altered based on the sensed data from Irvin’s sensor includes terminating audio output provided to a user.</p> <p>A POSITA would have been motivated to combine the teachings of the DHA system with Anderson and Abreu or Irvin to meet the claim limitations with a reasonable expectation of success for the reasons described above with respect to claim [1a] and [23a], and [35b].</p> <p>To the extent that the DHA system, alone or in view of Anderson and Abreu or Irvin, does not disclose any elements of this claim limitation, it would have been obvious to a person having ordinary skill in the art to combine, with the teachings of DHA systems, those elements, as known to them or as disclosed by other prior art identified in the cover pleading.</p>
<p>42[pre]. A personal audio system as recited in claim 24, wherein the head-worn wearable electronic device comprises:</p>	<p>DHA systems, alone or in view of Anderson and/or Enzmann, expressly or inherently discloses “[a] personal audio system as recited in claim 24, wherein the head-worn wearable electronic device comprises.”</p> <p><i>See</i> above at claim 1[pre].</p> <p>To the extent that the DHA system, alone or in view of Anderson, does not disclose any elements of this claim limitation, it would have been obvious to a person having ordinary skill in the art to combine, with the</p>

Claims	Prior Art
	<p>teachings of DHA systems, those elements, as known to them or as disclosed by other prior art identified in the cover pleading.</p>
<p>[42a] a first head-worn structure including at least a first speaker of the pair of speakers, the first speaker being wearable proximate to a first ear of the user; and</p>	<p>DHA systems, alone or in view of Anderson, expressly or inherently discloses “a first head-worn structure including at least a first speaker of the pair of speakers, the first speaker being wearable proximate to a first ear of the user.”</p> <p>Phonak – Perseo Brochure at p. 8</p> 

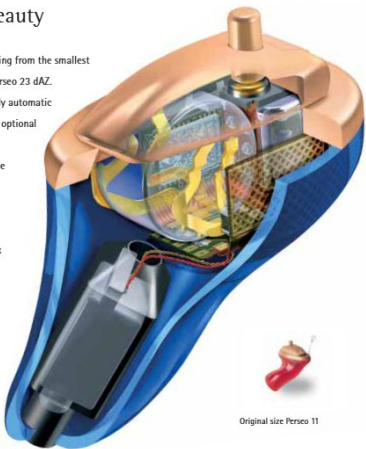
Phonak – Perseo Brochure at Page 22

The beauty of high technology

Perseo[®] ITEs – hidden beauty

The Perseo family includes a full line of ITEs ranging from the smallest Perseo 11 CIC to the elegant dual-microphone Perseo 23 dAZ. All CIC and ITE instruments, like the BTEs, are fully automatic and have **PersonalLogic**. A TacTronic switch is optional on all ITE models and easily added after fitting by exchanging battery covers. This is a very simple and practical means of providing extra control and additional programs.

In many countries, ITEs are available with Phonak eShell. These electronically created shells are hypoallergenic, exceptionally durable and more easily retained in the wearer's ear.



Original size Perseo 11

Valente - Effect of type of noise and loudspeaker array on the performance at 398

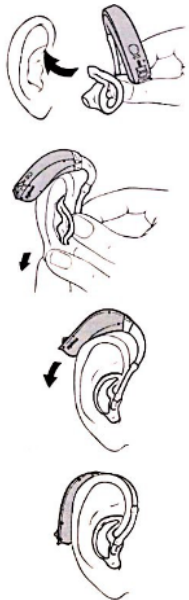
Differences in performance between omnidirectional and directional microphones were evaluated between two loudspeaker conditions (single loudspeaker at 180°; diffuse using eight loudspeakers set 45° apart) and two types of noise (steady-state HINT noise; R-Space[™] restaurant noise). Twenty-five participants were fit bilaterally with Phonak Perseo hearing aids using the manufacturer's recommended procedure. After wearing the hearing aids for one week, the

Claims

Prior Art

Phonak – Collections 2002 Product List at Page 8

	Claro 11 (CIC)	Claro 21 dAZ	Claro 22
ITE instruments			
			
Peak values G max./OSPL 90	45 dB/117 dB SPL	58 dB/117 dB SPL	58 dB/117 dB SPL
Features	World's first digital microphone TacTronic switch standard	In-the-Ear hearing computer with Adaptive digital Audio/oom (dAZ)	Semi-modular Claro In-the-Canal instrument
Description	Digital Perception Processing (DPP) AutoSelect (AS) Fine-scale Noise Canceler (FNC) 3 hearing programs + "mute" Battery size 10 TacTronic program selection switch	Digital Perception Processing (DPP) AutoSelect (AS) Fine-scale Noise Canceler (FNC) Adaptive digital Audio/oom (dAZ) 3 hearing programs + "mute" Battery size 312 Telecoil	Digital Perception Processing (DPP) AutoSelect (AS) Fine-scale Noise Canceler (FNC) 3 hearing programs + "mute" Battery size 312 Telecoil
Options	Rotary volume control potentiometer Colors: Pink, tan, cocoa, medium brown	Remote controls WatchPlot and SoundPlot FM TacTronic program selection switch On/Off switch (battery cover) Colors: Pink, tan, cocoa, medium brown	Remote controls WatchPlot and SoundPlot FM TacTronic program selection switch On/Off switch (battery cover) Colors: Pink, tan, cocoa, medium brown
Programming	Loudness Perception Profile (LPP) with Active Response Unit (ARU) Fitting formulas: Phonak Digital, Phonak Digital Ski Slope and DSL (i/o)	Loudness Perception Profile (LPP) with Active Response Unit (ARU) Fitting formulas: Phonak Digital, Phonak Digital Ski Slope and DSL (i/o)	Loudness Perception Profile (LPP) with Active Response Unit (ARU) Fitting formulas: Phonak Digital, Phonak Digital Ski Slope and DSL (i/o)

Claims	Prior Art
	<p data-bbox="583 175 1409 207">Oticon – DigiFocus Compact User Instructional Brochure at p. 6</p> <div data-bbox="632 224 1224 820">  <p data-bbox="850 329 1167 354">Putting on your hearing instruments</p> <ul data-bbox="850 358 1224 690" style="list-style-type: none"> • Switch off the hearing instrument. • Hold the plastic tube as shown. • Place the tip of the earmould in your ear canal, twisting slightly, making sure that the top part of the earmould (A) is pushed behind and under the fold (B) of your ear. • Gently pull down your ear lobe and press the earmould in the direction of the ear canal. • With the earmould positioned correctly in your ear, place the hearing instrument behind your ear by lifting its lower part and sliding it over the top of your ear. • Switch on the hearing instrument. </div> <p data-bbox="583 841 1077 873">Oticon – DigiFocus Description at p. 1</p> <p data-bbox="636 898 1312 966">DigiFocus is available both as behind-the-ear models (BTE) in colours which match the hair tone, and as in-the-ear (ITE) models which fit discreetly into the ear.</p> <div data-bbox="636 992 1320 1315">  <p data-bbox="716 1289 1241 1313"><i>DigiFocus Compact, Canal, Half Shell and Low Profile.</i></p> </div> <p data-bbox="583 1344 945 1377">Anderson '783 at 3:52–4:14</p> <p data-bbox="632 1382 1986 1450">One preferred embodiment of the invention herein described is shown in block diagram form in FIG. 1. In this figure, the earpiece 10 shown worn in the ear 11 uses a standard completely-in-the-canal (CIC)</p>

Claims	Prior Art
	<p>housing, well known in the art, although many other earpiece housing types known in the art (e.g., behind-the-ear, or BTE) can be used. The earpiece shown in FIG. 1 comprises a microphone 12, RF transceiver 13 with an antenna that doubles as an earpiece extractor 14, and speaker 15. Although the dual-use antenna/extractor arrangement shown in FIG. 1 is often desirable, a separate antenna or antennae can be used in conjunction with many different types of earpiece extractors as appropriate for a particular application. Antennae and extractors can be disguised as ear hair or jewelry (e.g., earrings), or conductive filaments can be permanently implanted in the ear cartilage to act as antennae using a process similar to ear piercing. The earpiece 10 communicates with a remote processor unit (RPU) 16 via a primary two-way RF link 17, although many other wireless link media (e.g., ultrasonic or infrared) may be used instead. The RPU is typically worn under clothing (e.g., carried in a pocket or purse), but may also be worn in plain sight (e.g., on a belt) if desired. The RPU 16 may be connected (via wired or wireless means 18) to optional secondary wireless link circuitry 19 that allows wireless communication between the RPU and other sources of information (e.g., the general subscriber telephone network) via a secondary wireless link. Note that the optional secondary wireless link circuitry 19 may or may not be contained within the RPU case 70 of FIG. 7.</p> <p>Anderson '783 at 4:30–33</p> <p>Although the system described here is monaural, it will be clear to those skilled in the art that a binaural (stereo) system can be constructed based on the same principles as the monaural system.</p> <p>Anderson '783 discloses that the earpiece of a hearing aid can be “worn in the ear,” as shown in FIG. 1, or behind the ear, or as part of headset 21, and that the earpiece includes a speaker. <i>See, e.g.</i>, Anderson at FIGS. 1, 2.</p> <p><i>See</i> above at claim [24a].</p> <p>To the extent that the DHA system, alone or in view of Anderson, does not disclose any elements of this claim limitation, it would have been obvious to a person having ordinary skill in the art to combine, with the teachings of DHA systems, those elements, as known to them or as disclosed by other prior art identified in the cover pleading.</p>
[42b] a second head-worn structure including at least a second speaker of the pair of	DHA systems, alone or in view of Anderson, expressly or inherently discloses “a second head-worn structure including at least a second speaker of the pair of speakers, the second speaker being wearable proximate to a second ear of the user.”

Claims	Prior Art
<p>speakers, the second speaker being wearable proximate to a second ear of the user.</p>	<p>Phonak – Perseo Brochure at p. 8</p>  <p>The image displays two black Phonak Perseo behind-the-ear hearing aids. Each device features a curved body with the 'PHONAK' logo and a clear, curved earhook. In the foreground, there are three custom-molded ear pieces: one blue, one red, and one tan. The tan pieces show the internal speaker and microphone components. The entire set is presented against a plain white background.</p>

Claims

Prior Art

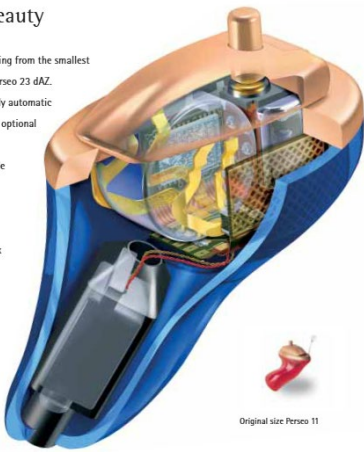
Phonak – Perseo Brochure at Page 22

The beauty of high technology

Perseo[®] ITEs – hidden beauty

The Perseo family includes a full line of ITEs ranging from the smallest Perseo 11 CIC to the elegant dual-microphone Perseo 23 dAZ. All CIC and ITE instruments, like the BTEs, are fully automatic and have **PersonalLogic**. A TacTronic switch is optional on all ITE models and easily added after fitting by exchanging battery covers. This is a very simple and practical means of providing extra control and additional programs.

In many countries, ITEs are available with Phonak eShell. These electronically created shells are hypoallergenic, exceptionally durable and more easily retained in the wearer's ear.

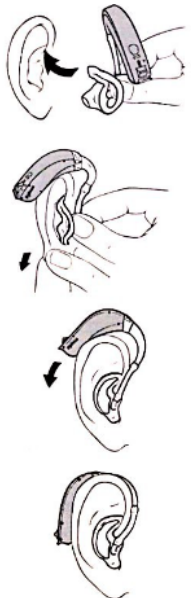


Original size Perseo 11

Valente - Effect of type of noise and loudspeaker array on the performance at 398

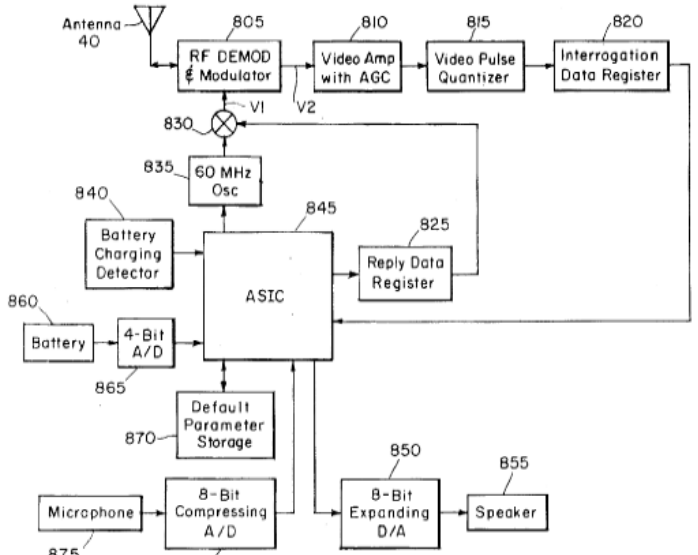
Differences in performance between omnidirectional and directional microphones were evaluated between two loudspeaker conditions (single loudspeaker at 180°; diffuse using eight loudspeakers set 45° apart) and two types of noise (steady-state HINT noise; R-Space[™] restaurant noise). Twenty-five participants were fit bilaterally with Phonak Perseo hearing aids using the manufacturer's recommended procedure. After wearing the hearing aids for one week, the

Claims	Prior Art																								
	<div>Phonak – Collections 2002 Product List at Page 8</div> <div><div><div>Claro 11 (CIC)</div><div>Claro 21 dAZ</div><div>Claro 22</div></div><div></div><div>ITE instruments</div><table><tr><td></td><td></td><td></td><td></td></tr><tr><td>Peak values G max./OSPL 90</td><td>45 dB/117 dB SPL</td><td>58 dB/117 dB SPL</td><td>58 dB/117 dB SPL</td></tr><tr><td>Features</td><td>World's first digital microphone TacTronic switch standard</td><td>In-the-Ear hearing computer with Adaptive digital Audio/oom (dAZ)</td><td>Semi-modular Claro In-the-Canal instrument</td></tr><tr><td>Description</td><td>Digital Perception Processing (DPP) AutoSelect (AS) Fine-scale Noise Canceler (FNC): 3 hearing programs + "mute" Battery size 10 TacTronic program selection switch</td><td>Digital Perception Processing (DPP) AutoSelect (AS) Fine-scale Noise Canceler (FNC): Adaptive digital Audio/oom (dAZ) 3 hearing programs + "mute" Battery size 312 Telecoil</td><td>Digital Perception Processing (DPP) AutoSelect (AS) Fine-scale Noise Canceler (FNC): 3 hearing programs + "mute" Battery size 312 Telecoil</td></tr><tr><td>Options</td><td>Rotary volume control potentiometer Colors: Pink, tan, cocoa, medium brown</td><td>Remote controls WatchPlot and SoundPlot FM TacTronic program selection switch On/off switch (battery cover) Colors: Pink, tan, cocoa, medium brown</td><td>Remote controls WatchPlot and SoundPlot FM TacTronic program selection switch On/off switch (battery cover) Colors: Pink, tan, cocoa, medium brown</td></tr><tr><td>Programming</td><td>Loudness Perception Profile (LPP) with Active Response Unit (ARU) Fitting formulas: Phonak Digital, Phonak Digital Ski Slope and DSL (i/o)</td><td>Loudness Perception Profile (LPP) with Active Response Unit (ARU) Fitting formulas: Phonak Digital, Phonak Digital Ski Slope and DSL (i/o)</td><td>Loudness Perception Profile (LPP) with Active Response Unit (ARU) Fitting formulas: Phonak Digital, Phonak Digital Ski Slope and DSL (i/o)</td></tr></table></div>					Peak values G max./OSPL 90	45 dB/117 dB SPL	58 dB/117 dB SPL	58 dB/117 dB SPL	Features	World's first digital microphone TacTronic switch standard	In-the-Ear hearing computer with Adaptive digital Audio/oom (dAZ)	Semi-modular Claro In-the-Canal instrument	Description	Digital Perception Processing (DPP) AutoSelect (AS) Fine-scale Noise Canceler (FNC): 3 hearing programs + "mute" Battery size 10 TacTronic program selection switch	Digital Perception Processing (DPP) AutoSelect (AS) Fine-scale Noise Canceler (FNC): Adaptive digital Audio/oom (dAZ) 3 hearing programs + "mute" Battery size 312 Telecoil	Digital Perception Processing (DPP) AutoSelect (AS) Fine-scale Noise Canceler (FNC): 3 hearing programs + "mute" Battery size 312 Telecoil	Options	Rotary volume control potentiometer Colors: Pink, tan, cocoa, medium brown	Remote controls WatchPlot and SoundPlot FM TacTronic program selection switch On/off switch (battery cover) Colors: Pink, tan, cocoa, medium brown	Remote controls WatchPlot and SoundPlot FM TacTronic program selection switch On/off switch (battery cover) Colors: Pink, tan, cocoa, medium brown	Programming	Loudness Perception Profile (LPP) with Active Response Unit (ARU) Fitting formulas: Phonak Digital, Phonak Digital Ski Slope and DSL (i/o)	Loudness Perception Profile (LPP) with Active Response Unit (ARU) Fitting formulas: Phonak Digital, Phonak Digital Ski Slope and DSL (i/o)	Loudness Perception Profile (LPP) with Active Response Unit (ARU) Fitting formulas: Phonak Digital, Phonak Digital Ski Slope and DSL (i/o)
Peak values G max./OSPL 90	45 dB/117 dB SPL	58 dB/117 dB SPL	58 dB/117 dB SPL																						
Features	World's first digital microphone TacTronic switch standard	In-the-Ear hearing computer with Adaptive digital Audio/oom (dAZ)	Semi-modular Claro In-the-Canal instrument																						
Description	Digital Perception Processing (DPP) AutoSelect (AS) Fine-scale Noise Canceler (FNC): 3 hearing programs + "mute" Battery size 10 TacTronic program selection switch	Digital Perception Processing (DPP) AutoSelect (AS) Fine-scale Noise Canceler (FNC): Adaptive digital Audio/oom (dAZ) 3 hearing programs + "mute" Battery size 312 Telecoil	Digital Perception Processing (DPP) AutoSelect (AS) Fine-scale Noise Canceler (FNC): 3 hearing programs + "mute" Battery size 312 Telecoil																						
Options	Rotary volume control potentiometer Colors: Pink, tan, cocoa, medium brown	Remote controls WatchPlot and SoundPlot FM TacTronic program selection switch On/off switch (battery cover) Colors: Pink, tan, cocoa, medium brown	Remote controls WatchPlot and SoundPlot FM TacTronic program selection switch On/off switch (battery cover) Colors: Pink, tan, cocoa, medium brown																						
Programming	Loudness Perception Profile (LPP) with Active Response Unit (ARU) Fitting formulas: Phonak Digital, Phonak Digital Ski Slope and DSL (i/o)	Loudness Perception Profile (LPP) with Active Response Unit (ARU) Fitting formulas: Phonak Digital, Phonak Digital Ski Slope and DSL (i/o)	Loudness Perception Profile (LPP) with Active Response Unit (ARU) Fitting formulas: Phonak Digital, Phonak Digital Ski Slope and DSL (i/o)																						

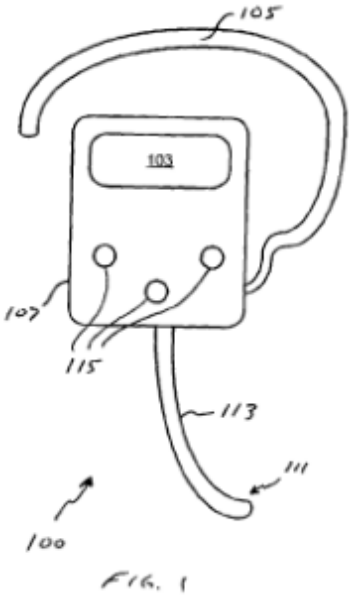
Claims	Prior Art
	<p data-bbox="583 175 1409 207">Oticon – DigiFocus Compact User Instructional Brochure at p. 6</p> <div data-bbox="632 224 1226 818">  <p data-bbox="848 329 1167 354">Putting on your hearing instruments</p> <ul data-bbox="848 358 1226 691" style="list-style-type: none"> • Switch off the hearing instrument. • Hold the plastic tube as shown. • Place the tip of the earmould in your ear canal, twisting slightly, making sure that the top part of the earmould (A) is pushed behind and under the fold (B) of your ear. • Gently pull down your ear lobe and press the earmould in the direction of the ear canal. • With the earmould positioned correctly in your ear, place the hearing instrument behind your ear by lifting its lower part and sliding it over the top of your ear. • Switch on the hearing instrument. </div> <p data-bbox="583 841 1077 873">Oticon – DigiFocus Description at p. 1</p> <p data-bbox="632 898 1314 967">DigiFocus is available both as behind-the-ear models (BTE) in colours which match the hair tone, and as in-the-ear (ITE) models which fit discreetly into the ear.</p>  <p data-bbox="709 1289 1241 1313"><i>DigiFocus Compact, Canal, Half Shell and Low Profile.</i></p> <p data-bbox="583 1344 947 1377">Anderson '783 at 3:52–4:14</p> <p data-bbox="632 1382 1986 1451">One preferred embodiment of the invention herein described is shown in block diagram form in FIG. 1. In this figure, the earpiece 10 shown worn in the ear 11 uses a standard completely-in-the-canal (CIC)</p>

Claims	Prior Art
	<p>housing, well known in the art, although many other earpiece housing types known in the art (e.g., behind-the-ear, or BTE) can be used. The earpiece shown in FIG. 1 comprises a microphone 12, RF transceiver 13 with an antenna that doubles as an earpiece extractor 14, and speaker 15. Although the dual-use antenna/extractor arrangement shown in FIG. 1 is often desirable, a separate antenna or antennae can be used in conjunction with many different types of earpiece extractors as appropriate for a particular application. Antennae and extractors can be disguised as ear hair or jewelry (e.g., earrings), or conductive filaments can be permanently implanted in the ear cartilage to act as antennae using a process similar to ear piercing. The earpiece 10 communicates with a remote processor unit (RPU) 16 via a primary two-way RF link 17, although many other wireless link media (e.g., ultrasonic or infrared) may be used instead. The RPU is typically worn under clothing (e.g., carried in a pocket or purse), but may also be worn in plain sight (e.g., on a belt) if desired. The RPU 16 may be connected (via wired or wireless means 18) to optional secondary wireless link circuitry 19 that allows wireless communication between the RPU and other sources of information (e.g., the general subscriber telephone network) via a secondary wireless link. Note that the optional secondary wireless link circuitry 19 may or may not be contained within the RPU case 70 of FIG. 7.</p> <p>Anderson '783 at 4:30–33</p> <p>Although the system described here is monaural, it will be clear to those skilled in the art that a binaural (stereo) system can be constructed based on the same principles as the monaural system.</p> <p>Anderson '783 discloses that the earpiece of a hearing aid can be “worn in the ear,” as shown in FIG. 1, or behind the ear, or as part of headset 21, and that the earpiece includes a speaker. <i>See, e.g.</i>, Anderson at FIGS. 1, 2. Anderson '783 discloses that a hearing aid system can be a “binaural (stereo) system.” <i>See id.</i> at 4:30–33. A POSITA would have understood from Anderson’s disclosures at 4:30–33 that Anderson’s hearing aid system can include a second earpiece, and thus a second speaker, to be worn in or near the user’s other ear, for example, when the earpiece is part of headset 21.</p> <p><i>See</i> above at claim [24a].</p> <p>To the extent that the DHA system, alone or in view of Anderson, does not disclose any elements of this claim limitation, it would have been obvious to a person having ordinary skill in the art to combine, with the teachings of DHA systems, those elements, as known to them or as disclosed by other prior art identified in the cover pleading.</p>

Claims	Prior Art
<p>43[pre]. A personal audio system as recited in claim 42, wherein the first head-worn structure comprises:</p>	<p>DHA systems, alone or in view of Anderson, expressly or inherently discloses “[a] personal audio system as recited in claim 42, wherein the first head-worn structure comprises.”</p> <p><i>See above at claims 1[pre]; [24a].</i></p> <p>A POSITA would have been motivated to combine the teachings of the DHA system with Anderson to meet the claim limitations with a reasonable expectation of success for the reasons described above with respect to claim [1a] and [26a].</p> <p>To the extent that the DHA system, alone or in view of Anderson, does not disclose any elements of this claim limitation, it would have been obvious to a person having ordinary skill in the art to combine, with the teachings of DHA systems, those elements, as known to them or as disclosed by other prior art identified in the cover pleading.</p>
<p>43[a] a first circuit board internal to the first head-worn structure,</p>	<p>DHA systems, alone or in view of Anderson and Jannard, expressly or inherently discloses “a first circuit board internal to the first head-worn structure.”</p> <p><i>See above at claim 31.</i></p> <p>A POSITA would have been motivated to combine the teachings of the DHA system with Anderson and Jannard to meet the claim limitations with a reasonable expectation of success for the reasons described above with respect to claim [1a] and [26a].</p> <p>To the extent that the DHA system, alone or in view of Anderson and Jannard, does not disclose any elements of this claim limitation, it would have been obvious to a person having ordinary skill in the art to combine, with the teachings of DHA systems, those elements, as known to them or as disclosed by other prior art identified in the cover pleading.</p>
<p>43[b] the first circuit board supporting at least a first integrated circuit and a second integrated circuit, the first integrated circuit and the second integrated circuit being physically attached to the first circuit board,</p>	<p>DHA systems, alone or in view of Anderson and Jannard, expressly or inherently discloses “the first circuit board supporting at least a first integrated circuit and a second integrated circuit, the first integrated circuit and the second integrated circuit being physically attached to the first circuit board.”</p> <p><i>See above at claim 31.</i></p> <p>A POSITA would have been motivated to combine the teachings of the DHA system with Anderson and Jannard to meet the claim limitations with a reasonable expectation of success for the reasons described above with respect to claim [1a] and [26a].</p> <p>To the extent that the DHA system, alone or in view of Anderson and Jannard, does not disclose any elements of this claim limitation, it would have been obvious to a person having ordinary skill in the art to combine,</p>

Claims	Prior Art
	<p>with the teachings of DHA systems, those elements, as known to them or as disclosed by other prior art identified in the cover pleading.</p>
<p>43[c] and the first integrated circuit being electrically connected to the second integrated circuit or another electronic component.</p>	<p>DHA systems, alone or in view of Anderson and Jannard, expressly or inherently discloses “and the first integrated circuit being electrically connected to the second integrated circuit or another electronic component.”</p> <p>Anderson '783 at FIG. 8</p>  <p>FIG. 8</p> <p>Anderson '783 at 3:40–41</p> <p>FIG. 8 is a block diagram showing details of an earpiece for a full-featured preferred embodiment of the invention.</p> <p>Anderson '783 at 4:30–33</p> <p>Although the system described here is monaural, it will be clear to those skilled in the art that a binaural (stereo) system can be constructed based on the same principles as the monaural system.</p> <p>Anderson '783 at 12:38–46</p> <p>Note that all earpiece functions (e.g., power control for the 60 MHz oscillator 835) are controlled by a</p>

Claims	Prior Art
	<p>custom application-specific integrated circuit (ASIC) 845. The ASIC 845 may be a low-power low-voltage device similar, e.g., to the ECI Semiconductor P576 custom programmable linear/digital gate array. Many of the functions that are shown as separate blocks for clarity in FIG. 8, e.g., reply data register 825, may be implemented within the ASIC 845.</p> <p>A POSITA would have been motivated to combine the teachings of the DHA system with Anderson and Jannard to meet the claim limitations with a reasonable expectation of success for the reasons described above with respect to claim [1a] and [26a].</p> <p>To the extent that the DHA system, alone or in view of Anderson and Jannard, does not disclose any elements of this claim limitation, it would have been obvious to a person having ordinary skill in the art to combine, with the teachings of DHA systems, those elements, as known to them or as disclosed by other prior art identified in the cover pleading.</p>
<p>44. A personal audio system as recited in claim 43, wherein the first circuit board is a flexible circuit board.</p>	<p>DHA systems, alone or in view of Anderson and Jannard, expressly or inherently discloses “[a] personal audio system as recited in claim 43, wherein the first circuit board is a flexible circuit board.”</p> <p><i>See above at claim 32.</i></p> <p>A POSITA would have been motivated to combine the teachings of the DHA system with Anderson and Jannard to meet the claim limitations with a reasonable expectation of success for the reasons described above with respect to claim [1a] and [26a].</p> <p>To the extent that the DHA system, alone or in view of Anderson and Jannard, does not disclose any elements of this claim limitation, it would have been obvious to a person having ordinary skill in the art to combine, with the teachings of DHA systems, those elements, as known to them or as disclosed by other prior art identified in the cover pleading.</p>
<p>45[pre]. A personal audio system as recited in claim 42, wherein the first head-worn structure comprises:</p>	<p>DHA systems, alone or in view of Anderson, expressly or inherently discloses “[a] personal audio system as recited in claim 42, wherein the first head-worn structure comprises.”</p> <p><i>See above at claims 1[pre]; [24a].</i></p> <p>A POSITA would have been motivated to combine the teachings of the DHA system with Anderson to meet the claim limitations with a reasonable expectation of success for the reasons described above with respect to claims [1a] and [1b].</p> <p>To the extent that the DHA system, alone or in view of Anderson, does not disclose any elements of this claim limitation, it would have been obvious to a person having ordinary skill in the art to combine, with the</p>

Claims	Prior Art
	<p>teachings of DHA systems, those elements, as known to them or as disclosed by other prior art identified in the cover pleading.</p>
<p>[45a] a user input device, the user input device being configured to receive a user input, the user input device is based on at least capacitive effect, wherein the user input being received induces at least a capacitive change to produce at least an electrical response.</p>	<p>DHA systems, alone or in view of Anderson, Bodley and Zadesky or Hollemans or Seshadri, expressly or inherently discloses “a user input device, the user input device being configured to receive a user input, the user input device is based on at least capacitive effect, wherein the user input being received induces at least a capacitive change to produce at least an electrical response.”</p> <p>Bodley '433 at FIG. 1</p>  <p>FIG. 1</p> <p>Bodley '433 at 2:7–8</p> <p>FIG. 1 is a frontal view of a generic headset that includes a display in accordance with the invention;</p> <p>Bodley '433 at 2:38–50</p> <p>Although headset 100 may use any of a variety of means of attaching to a user (e.g., headband, earpiece, etc.) and thus allowing hands-free operation of the headset, preferably headset 100 includes an earpiece 105 (e.g., an ear hook) that is coupled to the user's ear, or a portion thereof. More preferably, earpiece 105 fits around a portion of the user's ear. Earpiece 105 is coupled to body 107 of the headset, body 107 housing the various electronic components necessary to provide functionality to headset 100. Attached to</p>

Claims	Prior Art
	<p>the back surface of body 107 is an output transducer 109 (e.g., a speaker). Also coupled to body 107 is an input transducer 111 (e.g., a microphone), which, in this embodiment, is located at the end of a boom member 113.</p> <p>Bodley '433 at 2:65–3:22</p> <p>Display 103 can be used to provide the user with various types of information depending upon the design goals for the specific headset. The design goals are primarily based on the expected headset use, desired headset functionality, display size, display resolution, allowable headset power drain, and the processing capabilities of the headset on-board processor (also referred to as a controller or an interface controller). Examples of information that can be displayed on display 103 are headset battery level, headset signal level, call status (e.g., dialing, in use, etc.), caller identification, current time/date, and timer information such as elapsed time associated with a particular call. Display 103 can also be configured to flash or otherwise visually indicate an incoming call. Assuming that the information is provided to the headset processor, display 103 can also display the battery level of the base unit (e.g., cell phone, dongle, etc.) and the signal strength of the base unit (e.g., cell phone). Additionally, display 103 can be used during headset configuration (e.g., to adjust or set the ring mode, ring loudness, headset volume, display mode, voice dialing, time, date, etc.) or to provide additional headset functionality (e.g., phone lists, text messages, dialing using virtual numeric keypad or phone list dialing). Additionally, display 103 can be touch sensitive, thus providing a touch screen for data input (e.g., phone numbers into a phone list, dialing using a virtual numeric keypad).</p> <p>Zadesky '230 at 2:16–59</p> <p>With regards to touch pads, mice and track balls, a Cartesian coordinate system is used to monitor the position of the finger, mouse and ball, respectively, as they are moved. The Cartesian coordinate system is generally defined as a two dimensional coordinate system (x, y) in which the coordinates of a point (e.g., position of finger, mouse or ball) are its distances from two intersecting, often perpendicular straight lines, the distance from each being measured along a straight line parallel to each other. For example, the x, y positions of the mouse, ball and finger may be monitored. The x, y positions are then used to correspondingly locate and move the input pointer on the display screen.</p> <p>To elaborate further, touch pads generally include one or more sensors for detecting the proximity of the finger thereto. By way of example, the sensors may be based on resistive sensing, surface acoustic wave sensing, pressure sensing (e.g., strain gauge), optical sensing, capacitive sensing and the like. The sensors are generally dispersed about the touch pad with each sensor representing an x, y position. In most cases, the sensors are arranged in a grid of columns and rows. Distinct x and y position signals, which control</p>

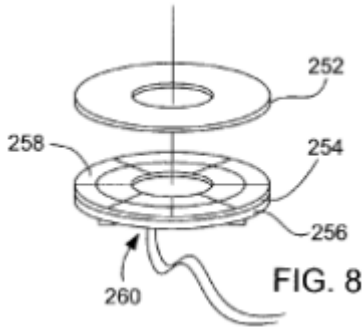
Claims	Prior Art
	<p>the x, y movement of a pointer device on the display screen, are thus generated when a finger is moved across the grid of sensors within the touch pad. For brevity sake, the remaining discussion will be held to the discussion of capacitive sensing technologies. It should be noted, however, that the other technologies have similar features.</p> <p>Capacitive sensing touch pads generally contain several layers of material. For example, the touch pad may include a protective shield, one or more electrode layers and a circuit board. The protective shield typically covers the electrode layer(s), and the electrode layer(s) is generally disposed on a front side of the circuit board. As is generally well known, the protective shield is the part of the touch pad that is touched by the user to implement cursor movements on a display screen. The electrode layer(s), on the other hand, is used to interpret the x, y position of the user's finger when the user's finger is resting or moving on the protective shield. The electrode layer (s) typically consists of a plurality of electrodes that are positioned in columns and rows so as to form a grid array. The columns and rows are generally based on the Cartesian coordinate system and thus the rows and columns correspond to the x and y directions.</p> <p>Zadesky '230 at 2:60–3:5</p> <p>The touch pad may also include sensing electronics for detecting signals associated with the electrodes. For example, the sensing electronics may be adapted to detect the change in capacitance at each of the electrodes as the finger passes over the grid. The sensing electronics are generally located on the backside of the circuit board. By way of example, the sensing electronics may include an application specific integrated circuit (ASIC) that is configured to measure the amount of capacitance in each of the electrodes and to compute the position of finger movement based on the capacitance in each of the electrodes. The ASIC may also be configured to report this information to the computing device.</p> <p>Zadesky '230 at 3:6–22</p> <p>Referring to FIG. 1, a touch pad 10 will be described in greater detail. The touch pad is generally a small rectangular area that includes a protective shield 12 and a plurality of electrodes 14 disposed underneath the protective shield layer 12. For case of discussion, a portion of the protective shield layer 12 has been removed to show the electrodes 14. Each of the electrodes 14 represents a different x, y position. In one configuration, as a finger 16 approaches the electrode grid 14, a tiny capacitance forms between the finger 16 and the electrodes 14 proximate the finger 16. The circuit board/sensing electronics measures capacitance and produces an x, y input signal 18 corresponding to the active electrodes 14 is sent to a host device 20 having a display screen 22. The x, y input signal 18 is used to control the movement of a cursor 24 on a display screen 22. As shown, the input pointer moves in a similar x, y direction as the detected x, y finger motion.</p>

Claims	Prior Art
	<p>Zadesky '230 at 12:31–41</p> <p>The capacitance in each of the electrodes 258 is measured by control circuitry 260 located on the backside of the circuit board 256. By detecting changes in capacitance at each of the electrodes 258, the control circuitry 260 can determine the angular location, direction, speed and acceleration of the finger as it is moved across the touch pad 250. The control circuitry 260 can also report this information in a form that can be used by a computing device. By way of example, the control circuitry may include an ASIC (application specific integrated circuit).</p> <p>Hollemans '490 at Abstract</p> <p>In a personal audio system (100), a remote controller (120) has a touch-sensitive area (122) and is worn in or by a human ear. A temporal pattern in the area (122), which is being touched, is detected and used for remotely controlling a device (110) for personal audio by means of a control signal (130). This prevents the hassle involved in finding, manipulating and operating a conventional remote control that is typically dangling somewhere along a wire.</p> <p>Hollemans '490 at 4:30–5:7</p> <p>The outer surface 121 is just the physical outer side of the controller 120. A part of this outer surface 121 is sensitive to being touched and is referred to as the touch-sensitive area 122. Optionally, the controller 120 comprises a second area that is sensitive to being touched, referred to as the further touch-sensitive area 123. The touch-sensitive area 122 can be realized in many ways, including a sensor for changes in capacity or induced voltages or currents, a grid of conducting parts between which changes in resistance are measured and detected, or an ordinary pressure sensor, button or temperature sensor. Another possibility is to measure and detect changes in either a passive or an active electromagnetic field. With the passive variant, the controller 120 relies on, for example, nearby power lines to generate measurable fields when being touched. With the active variant, the controller 120 generates an electromagnetic field of its own.</p> <p>Hollemans '490 at 8:19–22</p> <p>In a third way, the touch-sensitive area is used for picking up a gesture made on the area with e.g. a stylus or a finger. The remote control may be extended for this purpose with means for detecting a gesture, e.g. a laser beetle, a small camera, an IR movement sensor, or a cross-capacitance sensor.</p>

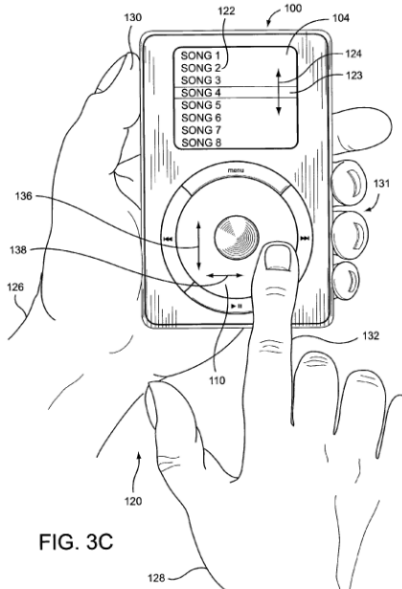
Claims	Prior Art
	<p>Hollemans '490 at 8:30–33</p> <p>The gesture detected may be a user defined gesture for activating a user defined function. The gesture detected may be used for user identification by matching with initials or an autograph. The gesture may extend into two dimensions, e.g. a planar gesture, or it may extend into three dimensions.</p> <p>Seshadri '839 at [0041]</p> <p>Wearable earpiece 12, as shown in FIG. 5, includes antenna 46, transmit/receive switch 122, receiver module 40, data recovery module 42, speaker module 44, transmitter module 120, input module 128 and display module 132. Receiver module 40, data recovery module 42 and speaker module 44 operate as discussed with reference to FIG. 4, with the exception of the data recovery module 42, to produce display information that is provided to display module 132. For instance, the received RF signal may include display information such as caller ID, command information, etc. which is separated by data recovery module 42 and provided to display module 132, which may be an LCD display, plasma display, etc.</p> <p>Seshadri '839 at [0042]</p> <p>Input module 128, which may be a keypad, touch screen, voice recognition circuit, or other like user interfaces, receives user commands and produces digital command messages 124 there from. Such digital command messages 124 includes, but are not limited to, packet size, synchronization information, frequency hopping initiation information, timeslot allocation information, link establishment information, piconet address information, fast-forward, play, pause, volume adjust, record, stop and rewind.</p> <p>A POSITA would have been motivated to combine the teachings of the DHA system with Anderson and Bodley or Seshadri or Hollemans or Zadesky to meet the claim limitations with a reasonable expectation of success for the reasons described above with respect to claims [1a] and [33]. In addition, a POSITA would have been motivated to combine the DHA system with the teachings of Bodley '433 or Seshadri to meet the claim limitation because Anderson '783 discloses that a user can control various functions and settings of an earpiece for a hearing aid system but faces the problem of how to adjust operation of the earpiece during use when the RPU, which is typically used to control operation of the earpiece, is disconnected. A POSITA would recognize Bodley '433 to provide a solution to this problem: providing the earpiece with a touch sensitive surface that can be used for controlling functions and settings of the earpiece. Anderson '783 and Bodley '433 are readily combinable with the DHA system because they share various structural and functional similarities, as explained in reference to claim 33. Moreover, a POSITA would be motivated to use Zadesky '230's capacitive sensing touch pad in the modified DHA system. Bodley '433 already discloses that a display for controlling settings and functions of an earpiece can include a touch sensitive surface. However,</p>

Claims	Prior Art
	<p>the specific details of the functionality of the touch sensitive surface are not disclosed in Bodley '433. A POSITA thus would have been motivated to look to other references to supply these missing details. A POSITA would have been motivated to look to Zadesky '230 or Hollemans or Seshadri at least because it is in the same general field of endeavor. As discussed above, each discloses head-worn devices for providing audio output to a user. Zadesky '230 discloses a media player that is “dedicated to processing media such as audio” and “allow[s] a user to listen to music.” Zadesky '230 and Seshadri are therefore from the same general field of endeavor as the DHA systems (such as the Phonak system), Anderson '783, and Bodley '433—providing audio output to a user. Thus, a POSITA would have been motivated to combine the teachings of these references with a reasonable expectation of success.</p> <p>To the extent that the DHA system, Anderson '783, Bodley '433, and Zadesky '230 or Hollemans or Seshadri does not disclose any elements of this claim limitation, it would have been obvious to a person having ordinary skill in the art to combine, with the teachings of DHA systems, those elements, as known to them or as disclosed by other prior art identified in the cover pleading.</p>
<p>46. A personal audio system as recited in claim 45, wherein the user input device is configured to distinguish at least a direction of the user input with respect to the user input device.</p>	<p>DHA systems, alone or in view of Anderson and Zadesky or Hollemans, expressly or inherently discloses “[a] personal audio system as recited in claim 45, wherein the user input device is configured to distinguish at least a direction of the user input with respect to the user input device.”</p> <p>Zadesky '230 at 2:16–59</p> <p>With regards to touch pads, mice and track balls, a Cartesian coordinate system is used to monitor the position of the finger, mouse and ball, respectively, as they are moved. The Cartesian coordinate system is generally defined as a two dimensional coordinate system (x, y) in which the coordinates of a point (e.g., position of finger, mouse or ball) are its distances from two intersecting, often perpendicular straight lines, the distance from each being measured along a straight line parallel to each other. For example, the x, y positions of the mouse, ball and finger may be monitored. The x, y positions are then used to correspondingly locate and move the input pointer on the display screen.</p> <p>To elaborate further, touch pads generally include one or more sensors for detecting the proximity of the finger thereto. By way of example, the sensors may be based on resistive sensing, surface acoustic wave sensing, pressure sensing (e.g., strain gauge), optical sensing, capacitive sensing and the like. The sensors are generally dispersed about the touch pad with each sensor representing an x, y position. In most cases, the sensors are arranged in a grid of columns and rows. Distinct x and y position signals, which control the x, y movement of a pointer device on the display screen, are thus generated when a finger is moved across the grid of sensors within the touch pad. For brevity sake, the remaining discussion will be held to the discussion of capacitive sensing technologies. It should be noted, however, that the other technologies have similar features.</p>

Claims	Prior Art
	<p>Capacitive sensing touch pads generally contain several layers of material. For example, the touch pad may include a protective shield, one or more electrode layers and a circuit board. The protective shield typically covers the electrode layer(s), and the electrode layer(s) is generally disposed on a front side of the circuit board. As is generally well known, the protective shield is the part of the touch pad that is touched by the user to implement cursor movements on a display screen. The electrode layer(s), on the other hand, is used to interpret the x, y position of the user's finger when the user's finger is resting or moving on the protective shield. The electrode layer (s) typically consists of a plurality of electrodes that are positioned in columns and rows so as to form a grid array. The columns and rows are generally based on the Cartesian coordinate system and thus the rows and columns correspond to the x and y directions.</p> <p>Zadesky '230 at 12:31–41</p> <p>The capacitance in each of the electrodes 258 is measured by control circuitry 260 located on the backside of the circuit board 256. By detecting changes in capacitance at each of the electrodes 258, the control circuitry 260 can determine the angular location, direction, speed and acceleration of the finger as it is moved across the touch pad 250. The control circuitry 260 can also report this information in a form that can be used by a computing device. By way of example, the control circuitry may include an ASIC (application specific integrated circuit).</p> <p>Hollemans '490 at 8:19–22</p> <p>In a third way, the touch-sensitive area is used for picking up a gesture made on the area with e.g. a stylus or a finger. The remote control may be extended for this purpose with means for detecting a gesture, e.g. a laser beetle, a small camera, an IR movement sensor, or a cross-capacitance sensor.</p> <p>Hollemans '490 at 8:30–33</p> <p>The gesture detected may be a user defined gesture for activating a user defined function. The gesture detected may be used for user identification by matching with initials or an autograph. The gesture may extend into two dimensions, e.g. a planar gesture, or it may extend into three dimensions.</p> <p>A POSITA would have been motivated to combine the teachings of the DHA system, Anderson '783, and Zadesky '230 or Hollemans '490 to meet the claim limitations with a reasonable expectation of success for the reasons described above with respect to claim 45.</p> <p>To the extent that the DHA system, Anderson '783 and Zadesky '230 or Hollemans '490 does not disclose any elements of this claim limitation, it would have been obvious to a person having ordinary skill in the art</p>

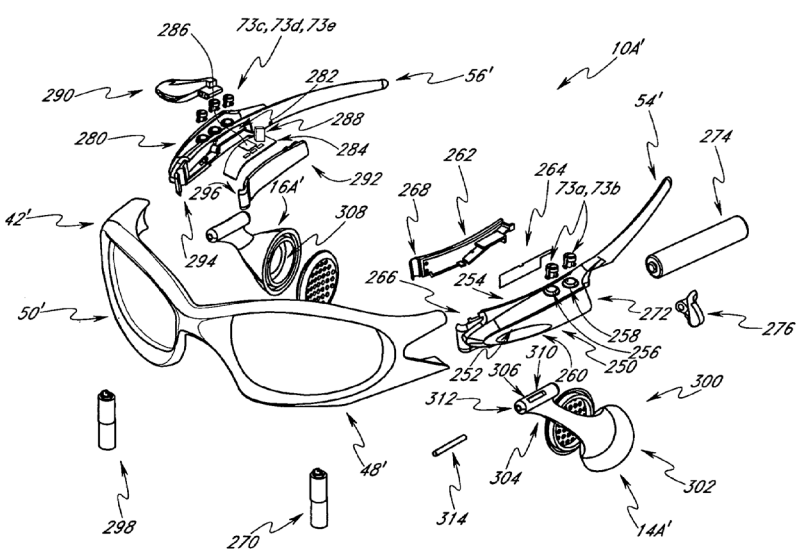
Claims	Prior Art
	<p>to combine, with the teachings of Anderson '783 and Zadesky '230, those elements, as known to them or as disclosed by other prior art identified in the cover pleading.</p>
<p>47. A personal audio system as recited in claim 46, wherein the user input device is configured to detect at least a sliding of a finger action with respect to the user input device.</p>	<p>DHA systems, alone or in view of Anderson and Zadesky or Hollemans, expressly or inherently discloses “[a] personal audio system as recited in claim 46, wherein the user input device is configured to detect at least a sliding of a finger action with respect to the user input device.”</p> <p>Zadesky '230 at FIG. 8</p>  <p>FIG. 8</p> <p>Zadesky '230 at 11:50–12:20</p> <p>FIG. 8 is a partially broken away perspective view of an annular capacitive touch pad 250, in accordance with one embodiment of the present invention. By way of example, the annular capacitive touch pad 250 may correspond to the touch pad of FIG. 2. The annular capacitive touch pad 250 is arranged to detect changes in capacitance as the user swirls an object such as a finger around the touch pad 250. The annular capacitive touch pad 250 is also arranged to detect changes in capacitance as the user moves their finger radially across the touch pad 250. The annular capacitive touch pad 250 is formed from various layers including at least a label layer 252, an electrode layer 254 and a circuit board 256. The label layer 252 is disposed over the electrode layer 254 and the electrode layer 254 is disposed over the circuit board 256. At least the label 252 and electrode layer 254 are annular such that they are defined by concentric circles, i.e., they have an inner perimeter and an outer perimeter. The circuit board 256 is generally a circular piece having an outer perimeter that coincides with the outer perimeter of the label 252 and electrode layer 254. It should be noted, however, that in some cases the circuit board 256 may be annular or the label 252 and electrode layer 254 may be circular.</p> <p>The label layer 252 serves to protect the underlayers and to provide a surface for allowing a finger to slide thereon. The surface is generally smooth so that the finger does not stick to it when moved. The label layer 252 also provides an insulating layer between the finger and the electrode layer 254. The electrode layer 254 includes a plurality of spatially distinct electrodes 258 that have positions based on the polar</p>

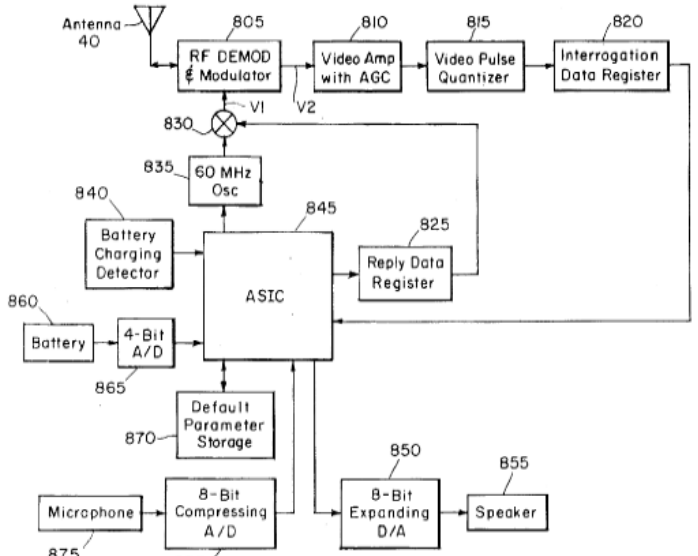
Claims	Prior Art
	<p>coordinate system. For instance, the electrodes 258 are positioned both angularly and radially on the circuit board 256 such that each of the electrodes 258 defines a distinct angular and radial position thereon. Any suitable number of electrodes 258 may be used. In most cases, it would be desirable to increase the number of electrodes 258 so as to provide higher resolution, i.e., more information can be used for things such as acceleration.</p> <p>Hollemans '490 at 8:19–22</p> <p>In a third way, the touch-sensitive area is used for picking up a gesture made on the area with e.g. a stylus or a finger. The remote control may be extended for this purpose with means for detecting a gesture, e.g. a laser beetle, a small camera, an IR movement sensor, or a cross-capacitance sensor.</p> <p>Hollemans '490 at 8:30–33</p> <p>The gesture detected may be a user defined gesture for activating a user defined function. The gesture detected may be used for user identification by matching with initials or an autograph. The gesture may extend into two dimensions, e.g. a planar gesture, or it may extend into three dimensions.</p> <p>A POSITA would have been motivated to combine the teachings of the DHA system with Anderson, Bodley, and Zadesky or Hollemans to meet the claim limitations with a reasonable expectation of success for the reasons described above with respect to claims [1a] and [45].</p> <p>To the extent that the DHA system, alone or in view of Anderson, Bodley, and Zadesky or Hollemans , does not disclose any elements of this claim limitation, it would have been obvious to a person having ordinary skill in the art to combine, with the teachings of DHA systems, those elements, as known to them or as disclosed by other prior art identified in the cover pleading</p>
<p>48. A personal audio system as recited in claim 45, wherein the user input device is configured to distinguish at least whether the user input pertains to a forward sliding action or a backward sliding action with respect to the user input device.</p>	<p>DHA systems, alone or in view of Anderson and Zadesky or Hollemans, expressly or inherently discloses “[a] personal audio system as recited in claim 45, wherein the user input device is configured to distinguish at least whether the user input pertains to a forward sliding action or a backward sliding action with respect to the user input device.”</p>

Claims	Prior Art
	<p>Zadesky '230 at FIG. 3C</p>  <p>FIG. 3C</p> <p>Zadesky '230 at 7:28–44</p> <p>The touch pad generally consists of a touchable outer surface 111 for receiving a finger for manipulation on the touch pad 110. Although not shown in FIG. 2, beneath the touchable outer surface 111 is a sensor arrangement. The sensor arrangement includes a plurality of sensors that are configured to activate as the finger passes over them. In the simplest case, an electrical signal is produced each time the finger passes a sensor. The number of signals in a given time frame may indicate location, direction, speed and acceleration of the finger on the touch pad, i.e., the more signals, the more the user moved his or her finger. In most cases, the signals are monitored by an electronic interface that converts the number, combination and frequency of the signals into location, direction, speed and acceleration information. This information may then be used by the media player 100 to perform the desired control function on the display screen 104.</p> <p>Zadesky '230 at 8:52–67</p> <p>FIGS. 3A-3C show the media player 100 of FIG. 2 being used by a user 120, in accordance with different embodiments of the invention. In all of these embodiments, the user 120 is linearly scrolling 104 (as shown by arrow 124) through a list of songs 122 displayed on the display screen via a slider bar 123. As</p>

Claims	Prior Art
	<p>shown, the media device 100 is comfortably held by one hand 126 while being comfortably addressed by the other hand 128. This configuration generally allows the user 120 to easily actuate the touch pad 110 with one or more fingers. For example, the thumb 130 and rightmost fingers 131 (or leftmost fingers if left handed) of the first hand 126 are used to grip the sides of the media player 100 while a finger 132 of the opposite hand 128 is used to actuate the touch pad 110. As shown, the entire top surface of the touch pad 110 is accessible to the user's finger 130.</p> <p>Zadesky '230 at 9:16–21</p> <p>Referring to FIG. 3C, and in accordance with one embodiment of the invention, the user 120 can slide his or her finger 132 substantially tangentially from all sides of the touch pad 110. For example, the touch pad 110 may be actuated forwards and backwards as shown by arrows 136 and side to side by arrows 138.</p> <p>Zadesky '230 at 11:50–12:20</p> <p>FIG. 8 is a partially broken away perspective view of an annular capacitive touch pad 250, in accordance with one embodiment of the present invention. By way of example, the annular capacitive touch pad 250 may correspond to the touch pad of FIG. 2. The annular capacitive touch pad 250 is arranged to detect changes in capacitance as the user swirls an object such as a finger around the touch pad 250. The annular capacitive touch pad 250 is also arranged to detect changes in capacitance as the user moves their finger radially across the touch pad 250. The annular capacitive touch pad 250 is formed from various layers including at least a label layer 252, an electrode layer 254 and a circuit board 256. The label layer 252 is disposed over the electrode layer 254 and the electrode layer 254 is disposed over the circuit board 256. At least the label 252 and electrode layer 254 are annular such that they are defined by concentric circles, i.e., they have an inner perimeter and an outer perimeter. The circuit board 256 is generally a circular piece having an outer perimeter that coincides with the outer perimeter of the label 252 and electrode layer 254. It should be noted, however, that in some cases the circuit board 256 may be annular or the label 252 and electrode layer 254 may be circular.</p> <p>The label layer 252 serves to protect the underlayers and to provide a surface for allowing a finger to slide thereon. The surface is generally smooth so that the finger does not stick to it when moved. The label layer 252 also provides an insulating layer between the finger and the electrode layer 254. The electrode layer 254 includes a plurality of spatially distinct electrodes 258 that have positions based on the polar coordinate system. For instance, the electrodes 258 are positioned both angularly and radially on the circuit board 256 such that each of the electrodes 258 defines a distinct angular and radial position thereon. Any suitable number of electrodes 258 may be used. In most cases, it would be desirable to</p>

Claims	Prior Art
	<p>increase the number of electrodes 258 so as to provide higher resolution, i.e., more information can be used for things such as acceleration.</p> <p>Hollemans '490 at 8:19–22 In a third way, the touch-sensitive area is used for picking up a gesture made on the area with e.g. a stylus or a finger. The remote control may be extended for this purpose with means for detecting a gesture, e.g. a laser beetle, a small camera, an IR movement sensor, or a cross-capacitance sensor.</p> <p>Hollemans '490 at 8:30–33 The gesture detected may be a user defined gesture for activating a user defined function. The gesture detected may be used for user identification by matching with initials or an autograph. The gesture may extend into two dimensions, e.g. a planar gesture, or it may extend into three dimensions.</p> <p>A POSITA would have been motivated to combine the teachings of the DHA system with Zadesky or Hollemans to meet the claim limitations with a reasonable expectation of success for the reasons described above with respect to claims [1a] and [45].</p> <p>To the extent that the DHA system, alone or in view of Zadesky or Hollemans, does not disclose any elements of this claim limitation, it would have been obvious to a person having ordinary skill in the art to combine, with the teachings of DHA systems, those elements, as known to them or as disclosed by other prior art identified in the cover pleading</p>
<p>[49a] A personal audio system as recited in claim 48, wherein the head-worn wearable electronic device comprises: a second circuit board internal to the second head-worn structure,</p>	<p>DHA systems, alone or in view of Jannard, expressly or inherently discloses “[a] personal audio system as recited in claim 48, wherein the head-worn wearable electronic device comprises: a second circuit board internal to the second head-worn structure.”</p>

Claims	Prior Art
	<p data-bbox="583 175 884 207">Jannard '526 at FIG. 3J</p>  <p data-bbox="919 782 1024 824"><i>FIG. 3J</i></p> <p data-bbox="583 844 894 876">Jannard '526 at 3:56–58</p> <p data-bbox="630 880 1990 948">FIG. 3J is a front, top, and left side perspective and exploded view of the wearable audio device illustrated in FIG. 3D.</p> <p data-bbox="583 1006 936 1039">Jannard '526 at 12:62–13:6</p> <p data-bbox="630 1042 1990 1295">A cover member 262 cooperates with the housing 250 to define the closed internal cavity 252. In the illustrated embodiment, the internal cavity 252 includes at least one compartment configured to receive an electronic circuit board 264 which includes at least one switch for each of the buttons 73a, 73b. In an exemplary but non-limiting embodiment, the board 264 can include two switches, one for each of the buttons 73a, 73b, which are configured to control a volume output from the speakers 14A', 16A'. The cover 262 can be attached to the ear stems 54A' with any type of fastener, such as, for example, but without limitation, screws, rivets, bolts, adhesive, and the like.</p> <p data-bbox="583 1354 911 1386">Jannard '526 at 13:50–55</p> <p data-bbox="630 1390 1990 1458">The internal cavity 282, in the illustrated embodiment, is configured to receive a printed circuit board 284. In the illustrated embodiment, the printed circuit board 284 includes one switch for each of the buttons</p>

Claims	Prior Art
	<p>73c, 73d, and 73e. Additionally, the printed circuit board 284 includes an audio file storage and playback device 286.</p> <p>A POSITA would have been motivated to combine the teachings of the DHA system with Jannard to meet the claim limitations with a reasonable expectation of success for the reasons described above with respect to claims [1a] and [26].</p> <p>To the extent that the DHA system, alone or in view of Jannard, does not disclose any elements of this claim limitation, it would have been obvious to a person having ordinary skill in the art to combine, with the teachings of DHA systems, those elements, as known to them or as disclosed by other prior art identified in the cover pleading.</p>
<p>[49b] the second circuit board supporting at least a third integrated circuit and a fourth integrated circuit, the third integrated circuit and the fourth integrated circuit being physically attached to the second circuit board, and</p>	<p>DHA systems, alone or in view of Anderson, expressly or inherently discloses “the second circuit board supporting at least a third integrated circuit and a fourth integrated circuit, the third integrated circuit and the fourth integrated circuit being physically attached to the second circuit board.”</p> <p>Anderson ’783 at FIG. 8</p>  <p>FIG. 8</p> <p>Anderson ’783 at 3:40–41</p> <p>FIG. 8 is a block diagram showing details of an earpiece for a full-featured preferred embodiment of the</p>

Claims	Prior Art
	<p>invention.</p> <p>Anderson '783 at 12:38–46</p> <p>Note that all earpiece functions (e.g., power control for the 60 MHz oscillator 835) are controlled by a custom application-specific integrated circuit (ASIC) 845. The ASIC 845 may be a low-power low-voltage device similar, e.g., to the ECI Semiconductor P576 custom programmable linear/digital gate array. Many of the functions that are shown as separate blocks for clarity in FIG. 8, e.g., reply data register 825, may be implemented within the ASIC 845.</p> <p>A POSITA would have been motivated to combine the teachings of the DHA system with Anderson to meet the claim limitations with a reasonable expectation of success for the reasons described above with respect to claims [1a].</p> <p>To the extent that the DHA system, alone or in view of Anderson, does not disclose any elements of this claim limitation, it would have been obvious to a person having ordinary skill in the art to combine, with the teachings of DHA systems, those elements, as known to them or as disclosed by other prior art identified in the cover pleading.</p>
<p>[49c] the third integrated circuit being electrically connected to the fourth integrated circuit or another electronic component.</p>	<p>DHA systems, alone or in view of Anderson and Jannard, expressly or inherently discloses “the third integrated circuit being electrically connected to the fourth integrated circuit or another electronic component.”</p> <p><i>See above at claim 43[c].</i></p> <p>A POSITA would have been motivated to combine the teachings of the DHA system with Anderson and Jannard to meet the claim limitations with a reasonable expectation of success for the reasons described above with respect to claims [1a] and [26].</p> <p>To the extent that the DHA system, alone or in view of Anderson and Jannard, does not disclose any elements of this claim limitation, it would have been obvious to a person having ordinary skill in the art to combine, with the teachings of DHA systems, those elements, as known to them or as disclosed by other prior art identified in the cover pleading.</p>
<p>50[pre]. A personal audio system as recited in claim 42, wherein the first head-worn structure of the</p>	<p>DHA systems, alone or in view of Anderson, expressly or inherently discloses “[a] personal audio system as recited in claim 42, wherein the first head-worn structure of the head-worn wearable electronic device comprises.”</p>

Claims	Prior Art
<p>head-worn wearable electronic device comprises:</p>	<p><i>See above at claim 1[pre]; 24.</i></p> <p>A POSITA would have been motivated to combine the teachings of the DHA system with Anderson to meet the claim limitations with a reasonable expectation of success for the reasons described above with respect to claims [1a] and 11.</p> <p>To the extent that the DHA system, alone or in view of Anderson and Nemirovski, does not disclose any elements of this claim limitation, it would have been obvious to a person having ordinary skill in the art to combine, with the teachings of DHA systems, those elements, as known to them or as disclosed by other prior art identified in the cover pleading.</p>
<p>[50a] a first microphone configured to pick up at least audio input from at least the user, the first microphone being at least partially in the first head-worn structure and is positioned to be directed towards the user when the first head-worn structure is being worn by the user; and</p>	<p>DHA systems, alone or in view of Anderson and Nemirovski, expressly or inherently discloses “a first microphone configured to pick up at least audio input from at least the user, the first microphone being at least partially in the first head-worn structure and is positioned to be directed towards the user when the first head-worn structure is being worn by the user.”</p> <p>Phonak – Perseo Brochure at Page 23</p>  <p>Perseo® BTEs – reliable beauty</p> <p>The Perseo BTE family includes a full choice of sizes, colors and power levels. All feature an elegant, award-winning design with a refined, ergonomic shape. They have a patented housing construction with bonded plastic and silicon to effectively keep out debris and moisture. Environmental protection extends further with patented, disposable microphone covers to ensure consistent, reliable performance, in all conditions. With dual-microphone systems this is particularly important since electronic matching cannot successfully compensate for performance reduction as a result of debris and moisture. A comprehensive range of colors for the housing also caters for specific tastes and trends.</p> <p>Original size New small Perseo BTE 111/211</p>

Claims	Prior Art
	<p data-bbox="583 175 1625 212">Valente - Effect of type of noise and loudspeaker array on the performance at 399</p> <p data-bbox="638 228 1020 493"> <i>Recently, Phonak introduced its digital hearing aid in which one feature is an adaptive directional microphone. This adaptive microphone reportedly ensures identical sensitivity and phase characteristics of the dual microphones and a mechanism to <i>automatically</i> change the polar pattern from omnidirectional to any directional pattern while compensating for the changes in low-frequency response associated with each pattern. This hearing aid can be programmed</i> </p> <p data-bbox="583 516 1512 553">Phonak – The Hearing Review Article re Perseo (Sep. 5, 2003) at Page 2</p> <p data-bbox="638 570 1224 862"> Perseo is a 20-channel digital hearing instrument. The default base-program is optimized for speech understanding in quiet environments. A second base-program with different frequency response settings and optional adaptive directional microphones provides additional support in noisy situations. The acoustic environment is continuously monitored and the appropriate base-program is automatically selected. This dual-base automatic function allows for independent fine tuning in one listening situation without compromising performance in the other. The automatic performance is designed to be easily personalized for the individual needs and preferences of the user. </p> <p data-bbox="583 885 1176 922">Phonak – Annual Report (2003–2004) at p. 12</p> <p data-bbox="638 938 1085 1328"> Within the FM business segment, the Group continued to leverage existing key technology into new applications. With SmartLink, Phonak has achieved another industry first by allowing hearing instrument users to access other communications systems. SmartLink combines a remote control, a transmitter, a directional microphone system and a Bluetooth link to mobile phone all in one single device. Phonak also presented EduLink, an innovative FM device to assist the learning development of children with Auditory Processing Disorder (APD). EduLink was perceived very positively within the medical community and recent classroom tests show favourable results. </p>

Claims

Prior Art

Phonak – Collections 2002 Product List at Page 7



Adaptive digital AudioZoom (dAZ)

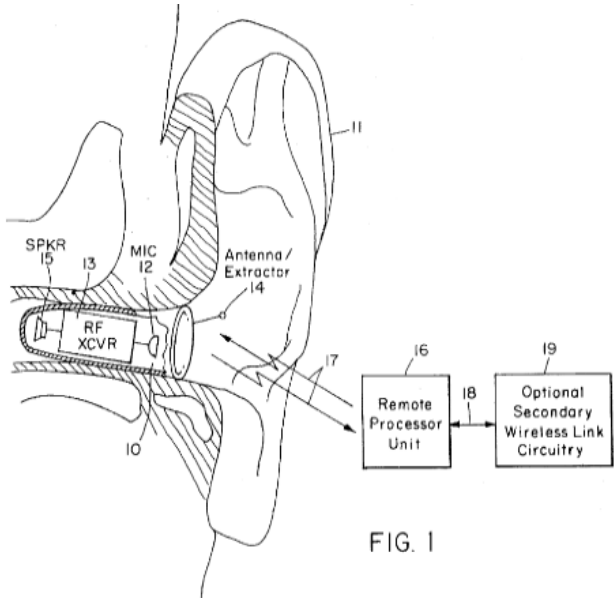


Fine-scale Noise Canceller (FNC)

<p>Active and adaptive separation of speech from noise The Speech-to-Noise Ratio (SNR) can only be optimised if the noise is spatially separated from the speech signal with directional microphone technology. As a pioneer in multi-microphone technologies, Phonak introduced AudioZoom, a feature which significantly improves listening and understanding in noisy situations, enhancing user satisfaction. Now, Phonak goes a step further in the refinement of multi-microphone technology with Adaptive digital AudioZoom: it seeks, finds, and suppresses noise automatically.</p>	<p>Fully automatic, fully adaptive With its unique spatial processing system, Phonak's Adaptive digital AudioZoom can detect noise not only from directly behind, but anywhere in the rear plane. Whether the hearing-impaired person is moving or the location of the noise source changes, the directional characteristic is constantly adapted to the listening situation. In this manner, the Adaptive digital AudioZoom automatically establishes the best possible SNR.</p>	<p>Twenty-band noise control Claro constantly monitors the Speech-to-Noise Ratio (SNR) in all 20 bands. If a poor SNR is measured in certain bands, the Fine-scale Noise Canceller will reduce gain only in those bands. The result of this precision noise control is a selective and intelligent improvement of the SNR with only a marginal influence on the speech signal.</p>	<p>Intelligent The magnitude of gain reduction is not constant. In frequency bands with higher SNRs and in bands which are important for speech intelligibility, gain is reduced only slightly. This intelligent control algorithm reduces the masking and distractive effects of background noise. Important speech information is preserved.</p>
	<p>Precise and selective All frequency bands are separately analysed and controlled. Gain is reduced selectively only in those bands where noise is present. This high-resolution analysis and gain control makes it possible to suppress noise with extremely high accuracy. Neighboring frequency bands without noise remain unaffected.</p>		

Lewis – Speech perception in noise: directional microphones (2004) at Page 429

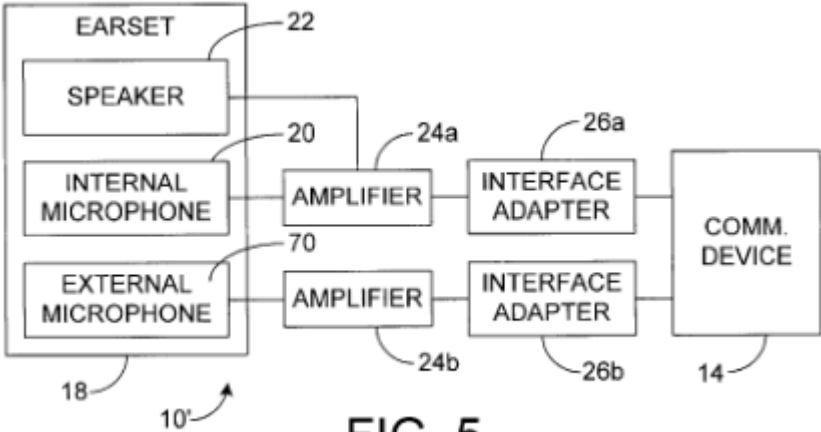
- (2) binaural digital Phonak Claro 311 dAZ BTE hearing aids alone in the omnidirectional microphone mode;
- (3) binaural digital Phonak Claro 311 dAZ BTE hearing aids alone in the directional (adaptive) microphone mode;
- (4) monaural digital Phonak Claro 311 dAZ BTE HA utilized with one Phonak Microlink FM receiver with the EM attenuated (FM-only mode) and one Phonak Claro 311 dAZ BTE HA in the omnidirectional microphone mode worn on the opposite ear; and
- (5) binaural digital Phonak Claro 311 dAZ BTE hearing aids utilized with binaural Phonak Microlink FM receivers with the EM attenuated (FM-only mode).

Claims	Prior Art
	<p data-bbox="583 175 894 207">Anderson '783 at FIG. 1</p>  <p data-bbox="995 748 1062 781">FIG. 1</p> <p data-bbox="583 846 945 878">Anderson '783 at 3:52–4:14</p> <p data-bbox="627 883 1986 1463">One preferred embodiment of the invention herein described is shown in block diagram form in FIG. 1. In this figure, the earpiece 10 shown worn in the ear 11 uses a standard completely-in-the-canal (CIC) housing, well known in the art, although many other earpiece housing types known in the art (e.g., behind-the-ear, or BTE) can be used. The earpiece shown in FIG. 1 comprises a microphone 12, RF transceiver 13 with an antenna that doubles as an earpiece extractor 14, and speaker 15. Although the dual-use antenna/extractor arrangement shown in FIG. 1 is often desirable, a separate antenna or antennae can be used in conjunction with many different types of earpiece extractors as appropriate for a particular application. Antennae and extractors can be disguised as ear hair or jewelry (e.g., earrings), or conductive filaments can be permanently implanted in the ear cartilage to act as antennae using a process similar to ear piercing. The earpiece 10 communicates with a remote processor unit (RPU) 16 via a primary two-way RF link 17, although many other wireless link media (e.g., ultrasonic or infrared) may be used instead. The RPU is typically worn under clothing (e.g., carried in a pocket or purse), but may also be worn in plain sight (e.g., on a belt) if desired. The RPU 16 may be connected (via wired or wireless means 18) to optional secondary wireless link circuitry 19 that allows wireless communication between the RPU and other sources of information (e.g., the general subscriber telephone network) via a secondary wireless link. Note that the optional secondary wireless link circuitry 19 may or may not be contained within the</p>

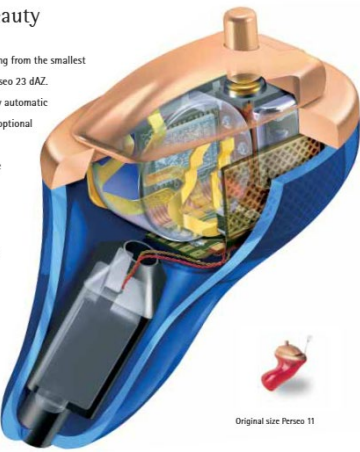
Claims	Prior Art
	<p>RPU case 70 of FIG. 7.</p> <p>Anderson '783 at 5:3–21</p> <p>During normal operation, speech signals from a nearby talker (and other signals in the ambient audio environment) are picked up by the earpiece wireless microphone 20 and transmitted to the RPU FM receiver 24. The RPU FM receiver 24 output level may be adjusted using the RPU FM receiver 24 volume control. The resulting electrical waveform representing signals from the nearby talker travels through the DPDT switch 25, which is set to the lower position as shown in FIG. 2, to the signal enhancer 26. The signal enhancer 26 may be, for example, a voice changer device that varies the pitch of a received speech signal according to the setting of pushbutton controls located on the RPU signal enhancer 26. The speech signal's pitch can then be raised or lowered as desired to help compensate for a user's hearing loss relative to a particular talker's voice characteristics. The modified speech signal travels from the RPU signal enhancer 26 to the RPU FM transmitter 27, and finally to the earpiece headset 21 receiver where the signal is converted to acoustic waves heard by the user.</p> <p>Anderson '783 at 6:8–25</p> <p>The user's voice (and other signals in the ambient audio environment) is picked up by the earpiece wireless microphone 20, transmitted to the RPU 23 FM receiver 24, sent through the DPDT switch 25 to the microphone input of the handset 28 (the handset's microphone transducer may be removed) and transmitted on the secondary wireless link 29 to, e.g., the general subscriber telephone network. Signals from the secondary wireless link 29 are received by the handset 28, sent from the speaker output of the handset 28 (the handset's speaker may be removed) via the DPDT switch 25 to the signal enhancer 26, FM transmitter 27 and finally to the earpiece headset 21 receiver where the signals are converted to acoustic waves heard by the user. Note that the RPU 23 may also include straight-forward connections to a "personal digital assistant" computer, or voice-operated devices connected to the FM receiver 24 output and not shown in FIG. 2.</p> <p>Anderson '783 at 22:11–61</p> <p>Stand-alone earpiece operation: If an earpiece cannot communicate with an RPU for any reason, stand-alone earpiece operation as a CIC hearing aid is generally desirable. Examples include occasions where the RPU must be turned off by pushing the RPU power-off key 78 (e.g., situations where operation of wireless transmitters is prohibited such as aircraft takeoffs and classified military briefings), RPU battery or other RPU electronics failure, misplacement of an RPU, and primary wireless link 17 interference. Note that a microphone located in the earpiece 10 is preferred to a body-worn microphone, since only an</p>

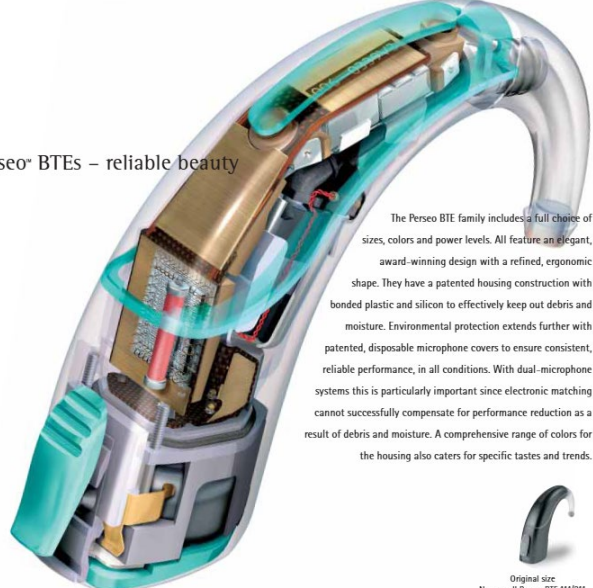
Claims	Prior Art
	<p>earpiece microphone 12 can be used to provide the desired default CIC hearing aid function (earpiece microphones also provide superior audio localization capabilities as the microphone moves with the user's head, and provide reduced wind and clothing noise compared with body-worn microphones). Parameters that control the stand-alone earpiece operation are stored in the earpiece CIC default parameter storage 870, which is accessed by the ASIC 845 to implement the desired stand-alone functions as described earlier. CIC default parameters include, but are not limited to, settings for the nominal and maximum allowable earpiece volume. Other information contained in the default parameter storage 870 includes, but is not limited to, the communication address and L/R (left/right) selection bit. To ensure that the CIC default parameter storage information cannot be changed inadvertently, changes are only allowed when the earpiece is not in use as, for example, when the earpiece battery 860 is being charged. A battery charging detector 840, that is sensitive to the high-level magnetic induction fields commonly used for earpiece battery charging, indicates to the ASIC 845 that the default parameter storage 870 values may be changed upon receipt of an appropriate command from the RPU. Specifically, to change information in the earpiece CIC default parameter storage 870, a user first activates the earpiece battery charging system, e.g., by placing the earpiece in a recharger/storage tray well known in the art. The RPU keyboard 946 is then used to enter a command to the DSP 948 that indicates the user's desire to change information in the default parameter storage 870, followed by the new information to be stored. The DSP 948 causes the RPU to transmit an interrogation with a reserved address code to which all earpieces being charged within the communication range of the RPU are receptive. The information that is to ultimately be loaded into the earpiece default parameter storage 870 is then transmitted by the RPU using the auxiliary bit field of many successive interrogations, said interrogations are received by the earpiece, and the desired information is loaded into the CIC default parameter storage 870 under control of the ASIC 845.</p> <p>Anderson '783 at 22:61–23:35</p> <p>Details of normal system operation: During normal operation of the hearing aid with wireless remote processor, sounds from the ambient environment are picked up by a microphone 12 in the earpiece 10 and transmitted with other information (e.g., data representing the earpiece battery 860 voltage as obtained by an A/D 865 converter and transmitted under control of the earpiece ASIC 845 in the auxiliary bit of many successive replies via the reply data register 825) over a primary two-way wireless link 17 to the RPU 16, where the audio signals are enhanced according to the user's needs by the RPU DSP 948. Signal enhancement is achieved via any of a number of techniques well known in the art that use the capabilities of the general-purpose RPU DSP 948. The enhanced audio signals may be combined with other information (e.g., a synthesized voice, generated within the DSP 948, warning that the aforementioned earpiece battery voltage is low) and transmitted from the RPU 16 over the primary wireless link 17 to the earpiece 10, where they are converted by a speaker 15 to sounds that can only be heard by the user 11. Art</p>

Claims	Prior Art
	<p>auxiliary speaker can be supplied in the RPU 16 as a peripheral 950 device that allows people other than the user to hear the same sounds as those produced by the earpiece speaker 15, if desired (as, e.g., when a person other than the user wishes to use the system's optional cellular telephone capabilities). A different, optional secondary two-way wireless link (shown as part of the RPU in FIG. 9 as secondary wireless link circuitry 944, or connected to the RPU by a wired or other, e.g. infrared, wireless link 18 and shown as optional secondary wireless link circuitry 19 in FIG. 1) can be used for communication between the RPU and a cellular telephone system or other sources of information. An RPU keyboard suitable for covert data entry 72, or voice recognition capabilities in the RPU (implemented, e.g., using any of a number of techniques well known in the art that use the general-purpose RPU DSP 948), can be used to control hearing aid parameters (such as amplification level set by the RPU DSP 948) and telephone dialing functions (e.g., automatic dialing of a predefined telephone number using information stored in the RPU DSP 948) in a manner that can be made imperceptible to a casual observer if desired.</p> <p>Anderson '783 at 23:36–42 Binaural, directional noise-canceling hearing aid with hearing protection: Two earpieces (one in each ear) and an RPU can be used in a binaural, directional noise-canceling hearing aid system with hearing protection features. Such a system provides some measure of protection for the residual hearing capability of a hearing-impaired user simultaneous with the hearing aid function.</p> <p>Anderson '783 at 25:29–30 The verbal commands are received as audio signals by the earpiece 10 and transmitted to the RPU 16.</p> <p>Anderson '783 at 26:54–65 Secondary link audio considerations: In some applications, a user may not wish to transmit signals from the ambient audio environment over the secondary wireless link. In such cases, an auxiliary earpiece microphone which favors the "bone conduction" path over the normal "air path" between the user's mouth and ear can be used. To save space, the auxiliary microphone may be housed in the same miniature case as the primary earpiece microphone. Alternatively, when two earpieces are used, unwanted ambient sources can be "nulled out" in the RPU DSP 948 by digital signal processing techniques well known in the art (e.g., the equivalent of speakerphone circuitry).</p> <p>Anderson '783 at 26:65–67 Note that such binaural signal processing techniques can also provide directivity and background noise</p>

Claims	Prior Art
	<p>cancellation for users</p> <p>Anderson '783 at 27:1–3 having total hearing loss in one ear and partial loss in the other ear, as discussed earlier in connection with binaural hearing aid system operation.</p> <p>Nemirovski '379 at FIG. 5</p>  <p style="text-align: center;">FIG. 5</p> <p>Nemirovski '379 at 2:59–60 FIG. 5 is a block diagram of an ear microphone assembly according to another aspect of the invention; and</p> <p>Nemirovski '379 at 3:6–23 The present invention is directed to a communication system, including a communications device and an ear microphone, or earset, assembly. The communications system allows a user to speak with a remotely located person. The ear microphone assembly, also referred to herein as an earset, includes a microphone and a speaker supported by a housing. The housing is retained by the ear of the user and allows for hands free use of the communication system while carrying on a conversation with the remotely located person. The microphone is directed towards the ear canal of the user and detects sounds emanating or coming out of the ear (or air pressure changes occurring within the ear) to accurately and reliably detect the speech of the user. The position of the microphone may also help to reduce detection of externally generated</p>

Claims	Prior Art
	<p>sounds. The ear microphone assembly allows for separation of a speech input device for a cellular, or wireless, telephone from an RF transceiver of the telephone.</p> <p>Nemirovski '379 at 9:45–60</p> <p>FIG. 5 is a block diagram that illustrates a system 10' for detecting speech of the user. The system 10' includes the earset 18 having the internal microphone 20, the speaker 22 and amplifier 24a as described in greater detail above. In addition, the earset 18 of the system 10' has the external microphone 70. As one skilled in the art will appreciate, the use of the external microphone 70 is optional. The external microphone 70 is disposed on the housing 30 (FIG. 3) and generally points away from the person wearing the earset 18. As a result, the external microphone 70 detects air pressure changes located outside the ear. These air pressure changes are mainly the result of externally generated noise (other people's speech, environmental noise, etc.) and, in part, to the speech of the person wearing the earset 18, when the speech travels from the mouth of the user and through the air to the external microphone 70.</p> <p>A POSITA would have been motivated to combine the teachings of the DHA system with Anderson and Nemirovski to meet the claim limitations with a reasonable expectation of success for the reasons described above with respect to claims [1a] and 11.</p> <p>To the extent that the DHA system, alone or in view of Anderson and Nemirovski, does not disclose any elements of this claim limitation, it would have been obvious to a person having ordinary skill in the art to combine, with the teachings of DHA systems, those elements, as known to them or as disclosed by other prior art identified in the cover pleading.</p>
<p>[50b] a second microphone configured to pick up at least ambient sound input associated with an environment in the vicinity of the first head-worn structure, the second microphone being positioned to be directed away from the user when the head-worn wearable electronic device is being worn by the user.</p>	<p>DHA systems, alone or in view of Anderson and Nemirovski, expressly or inherently discloses “a second microphone configured to pick up at least ambient sound input associated with an environment in the vicinity of the first head-worn structure, the second microphone being positioned to be directed away from the user when the head-worn wearable electronic device is being worn by the user.”</p>

Claims	Prior Art
	<p>Phonak – Perseo Brochure at Page 22</p> <h1 data-bbox="655 224 1255 376">The beauty of high technology</h1> <p data-bbox="703 376 961 402">Perseo[®] ITEs – hidden beauty</p> <div data-bbox="630 487 655 535">2</div> <p data-bbox="703 425 1003 587">The Perseo family includes a full line of ITEs ranging from the smallest Perseo 11 CIC to the elegant dual-microphone Perseo 23 dAZ. All CIC and ITE instruments, like the BTEs, are fully automatic and have PersonalLogic. A TacTronic switch is optional on all ITE models and easily added after fitting by exchanging battery covers. This is a very simple and practical means of providing extra control and additional programs.</p> <p data-bbox="703 613 919 695">In many countries, ITEs are available with Phonak eShell. These electronically created shells are hypoallergenic, exceptionally durable and more easily retained in the wearer's ear.</p>  <p data-bbox="1150 787 1234 799">Original size Perseo 11</p>

Claims	Prior Art
	<p data-bbox="583 175 1062 207">Phonak – Perseo Brochure at Page 23</p> <div data-bbox="625 224 1260 824">  <p data-bbox="636 383 905 407">Perseo® BTEs – reliable beauty</p> <p data-bbox="968 443 1255 695">The Perseo BTE family includes a full choice of sizes, colors and power levels. All feature an elegant, award-winning design with a refined, ergonomic shape. They have a patented housing construction with bonded plastic and silicon to effectively keep out debris and moisture. Environmental protection extends further with patented, disposable microphone covers to ensure consistent, reliable performance, in all conditions. With dual-microphone systems this is particularly important since electronic matching cannot successfully compensate for performance reduction as a result of debris and moisture. A comprehensive range of colors for the housing also caters for specific tastes and trends.</p> <p data-bbox="1140 800 1247 824">Original Size New small Perseo BTE 111/211</p> </div> <p data-bbox="583 846 1627 878">Valente - Effect of type of noise and loudspeaker array on the performance at 399</p> <p data-bbox="636 898 1022 1157">Recently, Phonak introduced its digital hearing aid in which one feature is an adaptive directional microphone. This adaptive microphone reportedly ensures identical sensitivity and phase characteristics of the dual microphones and a mechanism to <i>automatically</i> change the polar pattern from omnidirectional to any directional pattern while compensating for the changes in low-frequency response associated with each pattern. This hearing aid can be programmed</p>

Claims	Prior Art
	<p data-bbox="583 175 1514 212">Phonak – The Hearing Review Article re Perseo (Sep. 5, 2003) at Page 2</p> <p data-bbox="636 228 1226 521">Perseo is a 20-channel digital hearing instrument. The default base-program is optimized for speech understanding in quiet environments. A second base-program with different frequency response settings and optional adaptive directional microphones provides additional support in noisy situations. The acoustic environment is continuously monitored and the appropriate base-program is automatically selected. This dual-base automatic function allows for independent fine tuning in one listening situation without compromising performance in the other. The automatic performance is designed to be easily personalized for the individual needs and preferences of the user.</p> <p data-bbox="583 545 1176 583">Phonak – Annual Report (2003–2004) at p. 12</p> <p data-bbox="636 599 1087 989">Within the FM business segment, the Group continued to leverage existing key technology into new applications. With SmartLink, Phonak has achieved another industry first by allowing hearing instrument users to access other communications systems. SmartLink combines a remote control, a transmitter, a directional microphone system and a Bluetooth link to mobile phone all in one single device. Phonak also presented EduLink, an innovative FM device to assist the learning development of children with Auditory Processing Disorder (APD). EduLink was perceived very positively within the medical community and recent classroom tests show favourable results.</p>

Claims

Prior Art

Phonak – Collections 2002 Product List at Page 7



Adaptive digital AudioZoom (dAZ)



Fine-scale Noise Canceller (FNC)

Active and adaptive separation of speech from noise
The Speech-to-Noise Ratio (SNR) can only be optimised if the noise is spatially separated from the speech signal with directional microphone technology. As a pioneer in multi-microphone technologies, Phonak introduced AudioZoom, a feature which significantly improves listening and understanding in noisy situations, enhancing user satisfaction. Now, Phonak goes a step further in the refinement of multi-microphone technology with Adaptive digital AudioZoom: it seeks, finds, and suppresses noise automatically.

Fully automatic, fully adaptive
With its unique spatial processing system, Phonak's Adaptive digital AudioZoom can detect noise not only from directly behind, but anywhere in the near plane. Whether the hearing-impaired person is moving or the location of the noise source changes, the directional characteristic is constantly adapted to the listening situation. In this manner, the Adaptive digital AudioZoom automatically establishes the best possible SNR.

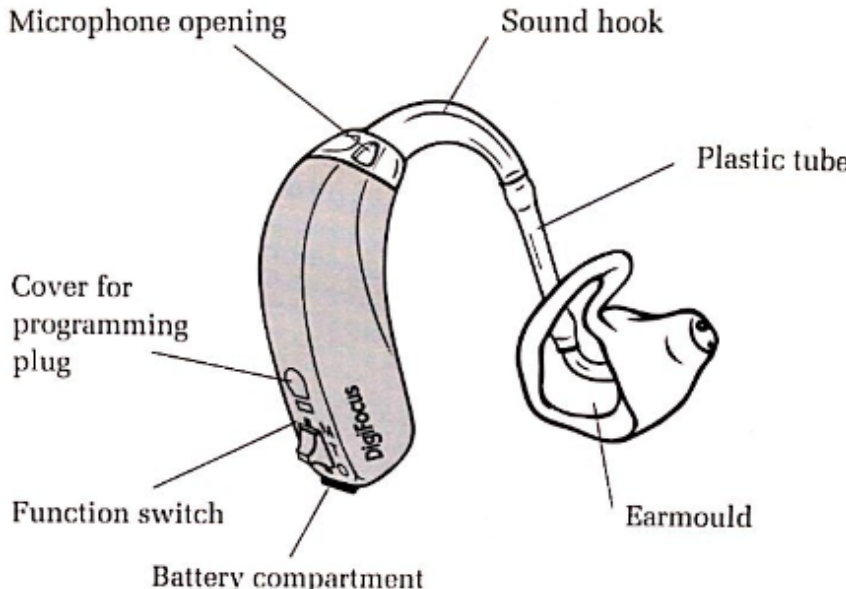
Twenty-band noise control
Claro constantly monitors the Speech-to-Noise Ratio (SNR) in all 20 bands. If a poor SNR is measured in certain bands, the Fine-scale Noise Canceller will reduce gain only in those bands. The result of this precision noise control is a selective and intelligent improvement of the SNR with only a marginal influence on the speech signal.

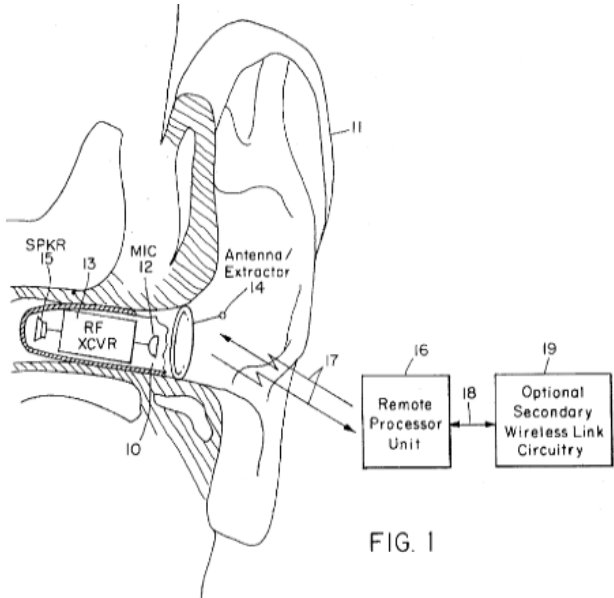
Intelligent
The magnitude of gain reduction is not constant. In frequency bands with higher SNRs and in bands which are important for speech intelligibility, gain is reduced only slightly. This intelligent control algorithm reduces the masking and distractive effects of background noise. Important speech information is preserved.

Precise and selective
All frequency bands are separately analysed and controlled. Gain is reduced selectively only in those bands where noise is present. This high-resolution analysis and gain control makes it possible to suppress noise with extremely high accuracy. Neighboring frequency bands without noise remain unaffected.

Lewis – Speech perception in noise: directional microphones (2004) at Page 429

- (2) binaural digital Phonak Claro 311 dAZ BTE hearing aids alone in the omnidirectional microphone mode;
- (3) binaural digital Phonak Claro 311 dAZ BTE hearing aids alone in the directional (adaptive) microphone mode;
- (4) monaural digital Phonak Claro 311 dAZ BTE HA utilized with one Phonak Microlink FM receiver with the EM attenuated (FM-only mode) and one Phonak Claro 311 dAZ BTE HA in the omnidirectional microphone mode worn on the opposite ear; and
- (5) binaural digital Phonak Claro 311 dAZ BTE hearing aids utilized with binaural Phonak Microlink FM receivers with the EM attenuated (FM-only mode).

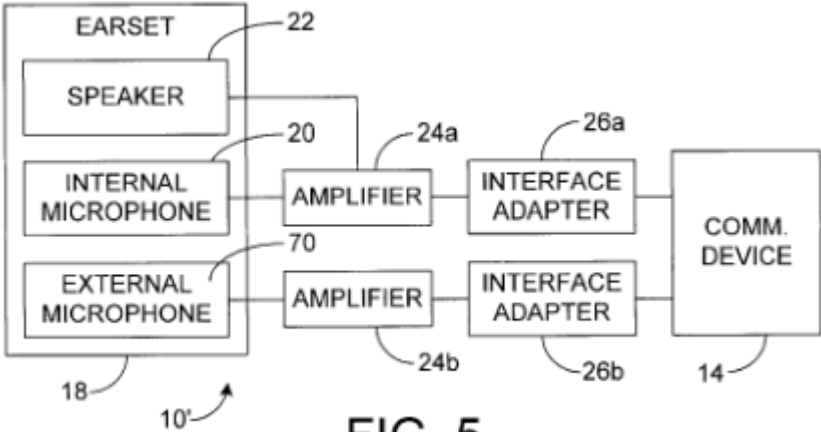
Claims	Prior Art
	<p data-bbox="583 175 1449 211">Oticon – DigiFocus Compact User Instructional Brochure at Page 4</p>  <p>The diagram illustrates the Oticon DigiFocus Compact hearing aid. It features a C-shaped body with a microphone opening at the top. A sound hook is attached to the upper part of the body. A plastic tube connects the body to an earmould. The body has a cover for a programming plug, a function switch, and a battery compartment. The earmould is shown fitting into the ear.</p>

Claims	Prior Art
	<p data-bbox="583 175 894 207">Anderson '783 at FIG. 1</p>  <p data-bbox="995 748 1062 781">FIG. 1</p> <p data-bbox="583 846 945 878">Anderson '783 at 3:52–4:14</p> <p data-bbox="627 883 1992 1463">One preferred embodiment of the invention herein described is shown in block diagram form in FIG. 1. In this figure, the earpiece 10 shown worn in the ear 11 uses a standard completely-in-the-canal (CIC) housing, well known in the art, although many other earpiece housing types known in the art (e.g., behind-the-ear, or BTE) can be used. The earpiece shown in FIG. 1 comprises a microphone 12, RF transceiver 13 with an antenna that doubles as an earpiece extractor 14, and speaker 15. Although the dual-use antenna/extractor arrangement shown in FIG. 1 is often desirable, a separate antenna or antennae can be used in conjunction with many different types of earpiece extractors as appropriate for a particular application. Antennae and extractors can be disguised as ear hair or jewelry (e.g., earrings), or conductive filaments can be permanently implanted in the ear cartilage to act as antennae using a process similar to ear piercing. The earpiece 10 communicates with a remote processor unit (RPU) 16 via a primary two-way RF link 17, although many other wireless link media (e.g., ultrasonic or infrared) may be used instead. The RPU is typically worn under clothing (e.g., carried in a pocket or purse), but may also be worn in plain sight (e.g., on a belt) if desired. The RPU 16 may be connected (via wired or wireless means 18) to optional secondary wireless link circuitry 19 that allows wireless communication between the RPU and other sources of information (e.g., the general subscriber telephone network) via a secondary wireless link. Note that the optional secondary wireless link circuitry 19 may or may not be contained within the</p>

Claims	Prior Art
	<p>RPU case 70 of FIG. 7.</p> <p>Anderson '783 at 5:3–21</p> <p>During normal operation, speech signals from a nearby talker (and other signals in the ambient audio environment) are picked up by the earpiece wireless microphone 20 and transmitted to the RPU FM receiver 24. The RPU FM receiver 24 output level may be adjusted using the RPU FM receiver 24 volume control. The resulting electrical waveform representing signals from the nearby talker travels through the DPDT switch 25, which is set to the lower position as shown in FIG. 2, to the signal enhancer 26. The signal enhancer 26 may be, for example, a voice changer device that varies the pitch of a received speech signal according to the setting of pushbutton controls located on the RPU signal enhancer 26. The speech signal's pitch can then be raised or lowered as desired to help compensate for a user's hearing loss relative to a particular talker's voice characteristics. The modified speech signal travels from the RPU signal enhancer 26 to the RPU FM transmitter 27, and finally to the earpiece headset 21 receiver where the signal is converted to acoustic waves heard by the user.</p> <p>Anderson '783 at 6:8–25</p> <p>The user's voice (and other signals in the ambient audio environment) is picked up by the earpiece wireless microphone 20, transmitted to the RPU 23 FM receiver 24, sent through the DPDT switch 25 to the microphone input of the handset 28 (the handset's microphone transducer may be removed) and transmitted on the secondary wireless link 29 to, e.g., the general subscriber telephone network. Signals from the secondary wireless link 29 are received by the handset 28, sent from the speaker output of the handset 28 (the handset's speaker may be removed) via the DPDT switch 25 to the signal enhancer 26, FM transmitter 27 and finally to the earpiece headset 21 receiver where the signals are converted to acoustic waves heard by the user. Note that the RPU 23 may also include straight-forward connections to a "personal digital assistant" computer, or voice-operated devices connected to the FM receiver 24 output and not shown in FIG. 2.</p> <p>Anderson '783 at 22:11–61</p> <p>Stand-alone earpiece operation: If an earpiece cannot communicate with an RPU for any reason, stand-alone earpiece operation as a CIC hearing aid is generally desirable. Examples include occasions where the RPU must be turned off by pushing the RPU power-off key 78 (e.g., situations where operation of wireless transmitters is prohibited such as aircraft takeoffs and classified military briefings), RPU battery or other RPU electronics failure, misplacement of an RPU, and primary wireless link 17 interference. Note that a microphone located in the earpiece 10 is preferred to a body-worn microphone, since only an</p>

Claims	Prior Art
	<p>earpiece microphone 12 can be used to provide the desired default CIC hearing aid function (earpiece microphones also provide superior audio localization capabilities as the microphone moves with the user's head, and provide reduced wind and clothing noise compared with body-worn microphones). Parameters that control the stand-alone earpiece operation are stored in the earpiece CIC default parameter storage 870, which is accessed by the ASIC 845 to implement the desired stand-alone functions as described earlier. CIC default parameters include, but are not limited to, settings for the nominal and maximum allowable earpiece volume. Other information contained in the default parameter storage 870 includes, but is not limited to, the communication address and L/R (left/right) selection bit. To ensure that the CIC default parameter storage information cannot be changed inadvertently, changes are only allowed when the earpiece is not in use as, for example, when the earpiece battery 860 is being charged. A battery charging detector 840, that is sensitive to the high-level magnetic induction fields commonly used for earpiece battery charging, indicates to the ASIC 845 that the default parameter storage 870 values may be changed upon receipt of an appropriate command from the RPU. Specifically, to change information in the earpiece CIC default parameter storage 870, a user first activates the earpiece battery charging system, e.g., by placing the earpiece in a recharger/storage tray well known in the art. The RPU keyboard 946 is then used to enter a command to the DSP 948 that indicates the user's desire to change information in the default parameter storage 870, followed by the new information to be stored. The DSP 948 causes the RPU to transmit an interrogation with a reserved address code to which all earpieces being charged within the communication range of the RPU are receptive. The information that is to ultimately be loaded into the earpiece default parameter storage 870 is then transmitted by the RPU using the auxiliary bit field of many successive interrogations, said interrogations are received by the earpiece, and the desired information is loaded into the CIC default parameter storage 870 under control of the ASIC 845.</p> <p>Anderson '783 at 22:61–23:35</p> <p>Details of normal system operation: During normal operation of the hearing aid with wireless remote processor, sounds from the ambient environment are picked up by a microphone 12 in the earpiece 10 and transmitted with other information (e.g., data representing the earpiece battery 860 voltage as obtained by an A/D 865 converter and transmitted under control of the earpiece ASIC 845 in the auxiliary bit of many successive replies via the reply data register 825) over a primary two-way wireless link 17 to the RPU 16, where the audio signals are enhanced according to the user's needs by the RPU DSP 948. Signal enhancement is achieved via any of a number of techniques well known in the art that use the capabilities of the general-purpose RPU DSP 948. The enhanced audio signals may be combined with other information (e.g., a synthesized voice, generated within the DSP 948, warning that the aforementioned earpiece battery voltage is low) and transmitted from the RPU 16 over the primary wireless link 17 to the earpiece 10, where they are converted by a speaker 15 to sounds that can only be heard by the user 11. Art</p>

Claims	Prior Art
	<p>auxiliary speaker can be supplied in the RPU 16 as a peripheral 950 device that allows people other than the user to hear the same sounds as those produced by the earpiece speaker 15, if desired (as, e.g., when a person other than the user wishes to use the system's optional cellular telephone capabilities). A different, optional secondary two-way wireless link (shown as part of the RPU in FIG. 9 as secondary wireless link circuitry 944, or connected to the RPU by a wired or other, e.g. infrared, wireless link 18 and shown as optional secondary wireless link circuitry 19 in FIG. 1) can be used for communication between the RPU and a cellular telephone system or other sources of information. An RPU keyboard suitable for covert data entry 72, or voice recognition capabilities in the RPU (implemented, e.g., using any of a number of techniques well known in the art that use the general-purpose RPU DSP 948), can be used to control hearing aid parameters (such as amplification level set by the RPU DSP 948) and telephone dialing functions (e.g., automatic dialing of a predefined telephone number using information stored in the RPU DSP 948) in a manner that can be made imperceptible to a casual observer if desired.</p> <p>Anderson '783 at 23:36–42 Binaural, directional noise-canceling hearing aid with hearing protection: Two earpieces (one in each ear) and an RPU can be used in a binaural, directional noise-canceling hearing aid system with hearing protection features. Such a system provides some measure of protection for the residual hearing capability of a hearing-impaired user simultaneous with the hearing aid function.</p> <p>Anderson '783 at 25:29–30 The verbal commands are received as audio signals by the earpiece 10 and transmitted to the RPU 16.</p> <p>Anderson '783 at 26:54–65 Secondary link audio considerations: In some applications, a user may not wish to transmit signals from the ambient audio environment over the secondary wireless link. In such cases, an auxiliary earpiece microphone which favors the "bone conduction" path over the normal "air path" between the user's mouth and ear can be used. To save space, the auxiliary microphone may be housed in the same miniature case as the primary earpiece microphone. Alternatively, when two earpieces are used, unwanted ambient sources can be "nulled out" in the RPU DSP 948 by digital signal processing techniques well known in the art (e.g., the equivalent of speakerphone circuitry).</p> <p>Anderson '783 at 26:65–67 Note that such binaural signal processing techniques can also provide directivity and background noise</p>

Claims	Prior Art
	<p>cancellation for users</p> <p>Anderson '783 at 27:1–3 having total hearing loss in one ear and partial loss in the other ear, as discussed earlier in connection with binaural hearing aid system operation.</p> <p>Nemirovski '379 at FIG. 5</p>  <p style="text-align: center;">FIG. 5</p> <p>Nemirovski '379 at 2:59–60 FIG. 5 is a block diagram of an ear microphone assembly according to another aspect of the invention; and</p> <p>Nemirovski '379 at 9:45–60 FIG. 5 is a block diagram that illustrates a system 10' for detecting speech of the user. The system 10' includes the earset 18 having the internal microphone 20, the speaker 22 and amplifier 24a as described in greater detail above. In addition, the earset 18 of the system 10' has the external microphone 70. As one skilled in the art will appreciate, the use of the external microphone 70 is optional. The external microphone 70 is disposed on the housing 30 (FIG. 3) and generally points away from the person wearing the earset 18. As a result, the external microphone 70 detects air pressure changes located outside the ear. These air pressure changes are mainly the result of externally generated noise (other people's speech, environmental noise, etc.) and, in part, to the speech of the person wearing the earset 18, when the speech</p>

Claims	Prior Art
	<p>travels from the mouth of the user and through the air to the external microphone 70.</p> <p>A POSITA would have been motivated to combine the teachings of the DHA system with Anderson and Nemirovski to meet the claim limitations with a reasonable expectation of success for the reasons described above with respect to claims [1a] and 11.</p> <p>To the extent that the DHA system, alone or in view of Anderson and Nemirovski, does not disclose any elements of this claim limitation, it would have been obvious to a person having ordinary skill in the art to combine, with the teachings of DHA systems, those elements, as known to them or as disclosed by other prior art identified in the cover pleading.</p>
<p>[51a] A personal audio system as recited in claim 42, wherein the first head-worn structure comprises: a first circuit board internal to the first head-worn structure, the first circuit board supporting at least a first integrated circuit, the first integrated circuit being physically attached to the first circuit board,</p>	<p>DHA systems, alone or in view of Anderson and Jannard, expressly or inherently discloses “[a] personal audio system as recited in claim 42, wherein the first head-worn structure comprises: a first circuit board internal to the first head-worn structure, the first circuit board supporting at least a first integrated circuit, the first integrated circuit being physically attached to the first circuit board.”</p> <p><i>See</i> above at claim 31.</p> <p>A POSITA would have been motivated to combine the teachings of the DHA system with Anderson and Jannard to meet the claim limitations with a reasonable expectation of success for the reasons described above with respect to claims [1a] and 71.</p> <p>To the extent that the DHA system, alone or in view of Anderson and Jannard does not disclose any elements of this claim limitation, it would have been obvious to a person having ordinary skill in the art to combine, with the teachings of DHA systems, those elements, as known to them or as disclosed by other prior art identified in the cover pleading.</p>
<p>[51b] wherein the first integrated circuit includes at least a processing unit, and</p>	<p>DHA systems, alone or in view of Anderson, expressly or inherently discloses “wherein the first integrated circuit includes at least a processing unit.”</p>

Claims

Prior Art

Phonak – Collections 2002 Product List at Page 6



Oticon – DigiFocus White Paper (1996) at Page 2

Manufacturers of hearing instruments have worked hard to overcome the troubles of size and battery consumption. With the thickness of digital microchips at only 1/100 of a human hair, this has now been accomplished in the world's first fully digital, at the ear hearing instrument - DigiFocus. Inside DigiFocus there is computing power available which equals that of a PC.

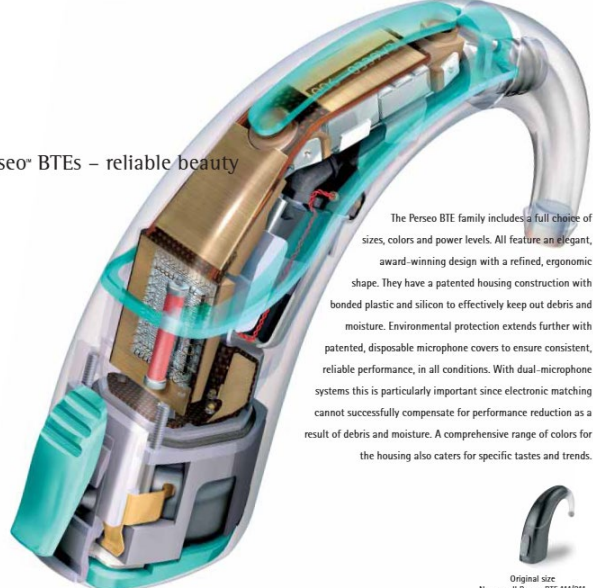
Anderson '783 at 12:38–46

Note that all earpiece functions (e.g., power control for the 60 MHz oscillator 835) are controlled by a custom application-specific integrated circuit (ASIC) 845. The ASIC 845 may be a low-power low-voltage device similar, e.g., to the ECI Semiconductor P576 custom programmable linear/digital gate array. Many of the functions that are shown as separate blocks for clarity in FIG. 8, e.g., reply data register 825, may be implemented within the ASIC 845.

Anderson '783 at 21:50–62

In the simplest case, the ASIC 845 reads data from the A/D 880, performs a multiplication of the data by a constant (read from the default parameter storage 870) to provide amplification, compares the result with maximum allowable values (read from the default parameter storage 870) and writes the result to the D/A

Claims	Prior Art
	<p>850. The ASIC 845 may perform other processing (e.g., adaptive volume control) as well. Although none of the remote processor and communication features are available in stand-alone mode, which uses pre-selected hearing aid parameter settings from the default parameter storage 870 in the earpiece, the earpiece provides a useful amplified audio path from the microphone to the speaker.</p> <p>A POSITA would have been motivated to combine the teachings of the DHA system with Anderson to meet the claim limitations with a reasonable expectation of success for the reasons described above with respect to claims [1a] and 71.</p> <p>To the extent that the DHA system, alone or in view of Anderson, does not disclose any elements of this claim limitation, it would have been obvious to a person having ordinary skill in the art to combine, with the teachings of DHA systems, those elements, as known to them or as disclosed by other prior art identified in the cover pleading.</p>
<p>[51c] wherein the processing unit is configured to: determine an ambient sound level; and</p>	<p>DHA systems, alone or in view of Jannard and Victorian, expressly or inherently discloses “wherein the processing unit is configured to: determine an ambient sound level.”</p> <p>Phonak – Perseo Brochure at Page 13</p> <div data-bbox="625 850 1484 896"> <h3>Enhanced sound quality and Speech-to-Noise Ratio</h3> </div> <div data-bbox="718 945 1064 971"> <p>Clarity Component Enhancement – DPP²</p> </div> <div data-bbox="718 982 1726 1203"> <p>Speech clarity depends on both its acoustic properties and the listener's auditory perception. Conventional spectral enhancement systems consider only the acoustics of sound, not how it will be perceived by the listener. Perseo goes well beyond this with DPP², the new generation of Digital Perception Processing. DPP² with Clarity Component Enhancement is individualized for each wearer. It analyzes the acoustic and the psychoacoustic characteristics of incoming sound in 20 critical bands to identify accurately the components that contribute to clarity. It then selectively amplifies them to optimize clarity for the wearer.</p> </div> 

Claims	Prior Art
	<p data-bbox="583 175 1062 207">Phonak – Perseo Brochure at Page 23</p> <div data-bbox="625 224 1260 824">  <p data-bbox="636 383 905 407">Perseo® BTEs – reliable beauty</p> <p data-bbox="968 443 1255 695">The Perseo BTE family includes a full choice of sizes, colors and power levels. All feature an elegant, award-winning design with a refined, ergonomic shape. They have a patented housing construction with bonded plastic and silicon to effectively keep out debris and moisture. Environmental protection extends further with patented, disposable microphone covers to ensure consistent, reliable performance, in all conditions. With dual-microphone systems this is particularly important since electronic matching cannot successfully compensate for performance reduction as a result of debris and moisture. A comprehensive range of colors for the housing also caters for specific tastes and trends.</p> <p data-bbox="1140 800 1247 824">Original size New small Perseo BTE 111/211</p> </div> <p data-bbox="583 846 1627 878">Valente - Effect of type of noise and loudspeaker array on the performance at 399</p> <p data-bbox="636 898 1022 1157">Recently, Phonak introduced its digital hearing aid in which one feature is an adaptive directional microphone. This adaptive microphone reportedly ensures identical sensitivity and phase characteristics of the dual microphones and a mechanism to <i>automatically</i> change the polar pattern from omnidirectional to any directional pattern while compensating for the changes in low-frequency response associated with each pattern. This hearing aid can be programmed</p>

Claims

Prior Art

Phonak – The Hearing Review Article re Perseo (Sep. 5, 2003) at Page 2

Perseo is a 20-channel digital hearing instrument. The default base-program is optimized for speech understanding in quiet environments. A second base-program with different frequency response settings and optional adaptive directional microphones provides additional support in noisy situations. The acoustic environment is continuously monitored and the appropriate base-program is automatically selected. This dual-base automatic function allows for independent fine tuning in one listening situation without compromising performance in the other. The automatic performance is designed to be easily personalized for the individual needs and preferences of the user.

Phonak – Collections 2002 Product List at Page 6

Elite


Digital Perception Processing (DPP)


AutoSelect (AS)

<p>New Phonak sound: unmatched clarity Digital Perception Processing (DPP) is based on a cochlear sound analysis across 20 critical bands – exactly like the analysis performed by the healthy human ear. But Clari gives a step further than previous overlapping multi-band filter systems. DPP models the active, load-independent cochlear signal processing in all 20 critical bands.</p>	<p>20 critical bands: overlapping and dynamic The cochlear filters are overlapping, but the degree of overlap is dependent on the frequency level, and bandwidth of the sound signal. This is precisely how Clari's filters work. Digital Perception Processing computes the cochlear excitation patterns and thus determines sound perception. DPP is the first to mimic the healthy human ear in this way.</p>	<p>Automatic, intelligent program selection Clari listens and automatically activates the optimal DPP algorithm for the acoustic environment. Program 1 with QuietAdapt for quiet situations or Program 2 with NoiseAdapt and noise control for noisy environments.</p> <p>In quiet situations QuietAdapt DPP is designed for clear hearing in quiet surroundings. QuietAdapt controls gain in 20 critical bands to correct for loudness and provide a natural sound quality. By preserving the perception patterns of complex sounds, QuietAdapt assures that low-level signal components are not over-</p>	<p>amplified to the point where they would adversely affect speech intelligibility. This allows Clari to deliver only those sound components that would be naturally audible.</p> <p>In noise NoiseAdapt DPP is the ideal program for maximizing audibility in noisy environments. The special NoiseAdapt algorithm works with the Fine-scale Noise Canceller – and if available with the Adaptive digital AudioZoom feature – to optimize speech understanding in noise.</p>
---	---	--	--



Claims

Prior Art

Phonak – Collections 2002 Product List at Page 7



Adaptive digital AudioZoom (dAZ)



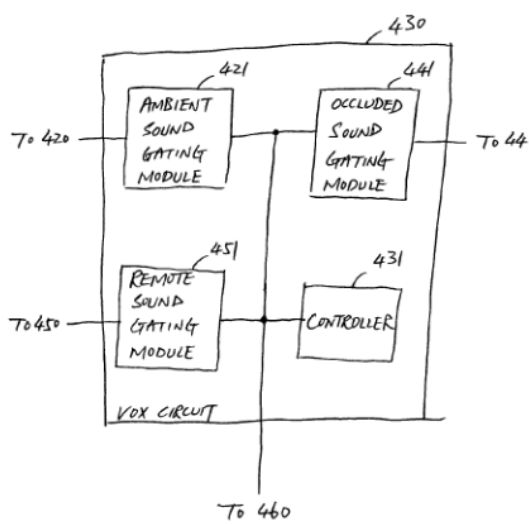
Fine-scale Noise Canceler (FNC)

<p>Active and adaptive separation of speech from noise The Speech-to-Noise Ratio (SNR) can only be optimised if the noise is spatially separated from the speech signal with directional microphone technology. As a pioneer in multi-microphone technologies, Phonak introduced AudioZoom, a feature which significantly improves listening and understanding in noisy situations, enhancing user satisfaction. Now, Phonak goes a step further in the refinement of multi-microphone technology with Adaptive digital AudioZoom: it seeks, finds, and suppresses noise automatically.</p>	<p>Fully automatic, fully adaptive With its unique spatial processing system, Phonak's Adaptive digital AudioZoom can detect noise not only from directly behind, but anywhere in the near plane. Whether the hearing-impaired person is moving or the location of the noise source changes, the directional characteristic is constantly adapted to the listening situation. In this manner, the Adaptive digital AudioZoom automatically establishes the best possible SNR.</p>	<p>Twenty-band noise control Claro constantly monitors the Speech-to-Noise Ratio (SNR) in all 20 bands. If a poor SNR is measured in certain bands, the Fine-scale Noise Canceler will reduce gain only in those bands. The result of this precision noise control is a selective and intelligent improvement of the SNR with only a marginal influence on the speech signal.</p>	<p>Intelligent The magnitude of gain reduction is not constant. In frequency bands with higher SNRs and in bands which are important for speech intelligibility, gain is reduced only slightly. This intelligent control algorithm reduces the masking and distractive effects of background noise. Important speech information is preserved.</p>
	<p>Precise and selective All frequency bands are separately analysed and controlled. Gain is reduced selectively only in those bands where noise is present. This high-resolution analysis and gain control makes it possible to suppress noise with extremely high accuracy. Neighboring frequency bands without noise remain unaffected.</p>		

Oticon – LA Time Article re THE GOODS: The Sound Wave of the Future (Nov. 10 1995) at 1
DigiFocus’ computer chips can remember different sound environments and adjust for them, depending on what listening situation the person is experiencing.

Jannard ’526 at 25:53–26:20

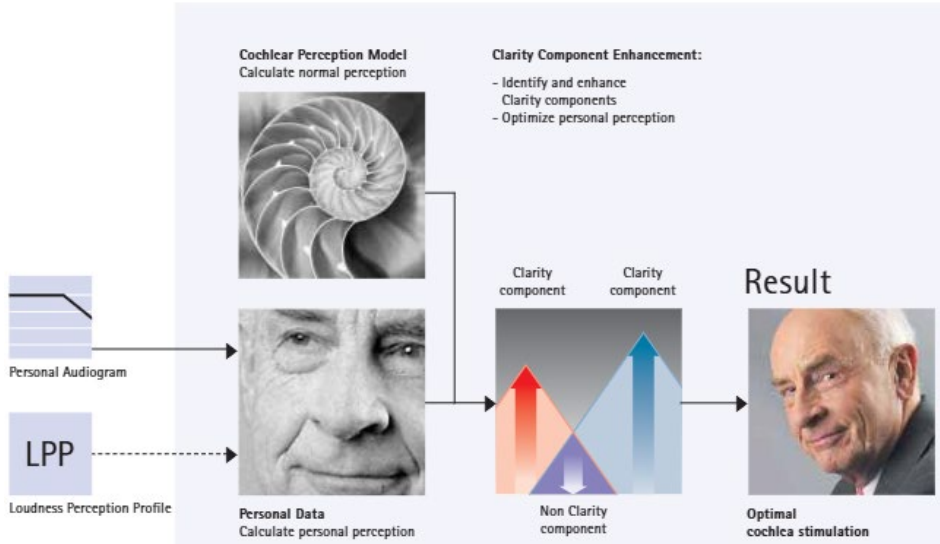
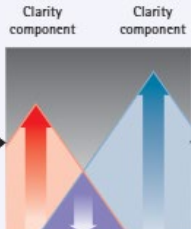
In another embodiment, a further modification of the audio devices 10, 10A, 10B, 10C, and 10D is provided wherein the audio devices include at least two banks of microphones, with one bank acting as a speaker of received and one bank providing an ambient noise-cancellation function. The microphone banks can be positioned at any suitable location or combination of locations (e.g., on the audio device, within the audio device, opposing sides of the audio device, or the like). In one embodiment, automatic switching of the speaking-microphone and noise-canceling-microphone banks' functions advantageously enhances ease of use. In a further embodiment, the microphone banks can be arranged in an array to be used in conjunction with algorithms to discern, reduce, and/or eliminate noise for the purpose of voice recognition. For example, in one embodiment, such microphone banks can include ASIC-based noise-canceling technology, such as is available in chips from Andrea Electronics Corporation (AEC), to enable voice recognition in ambient noise up to about 130 Db or more. In another embodiment, microphone

Claims	Prior Art
	<p>banks can be arranged in any suitable combination of linear or non-linear arrays to be used in conjunction with algorithms to discern, reduce, and/or eliminate noise in any suitable manner. In another embodiment, audio/proximity sensors can advantageously trigger the appropriate functionality in a specific bank. In another embodiment, a noise-canceling microphone can be provided in connection with a cord or other microphones described above. For example, without limitation, a series of miniature microphones can be supported down a cord from the audio device, separated by desired distances, and aimed in different directions. In another implementation, one or more of the microphones can be for verbal input from the user, and one or more others of the microphones, or the same microphone, can also be for noise-cancellation purposes.</p> <p>Victorian '313 at FIG. 4B</p>  <p>FIG. 4B</p> <p>Victorian '313 at [0021] FIG. 4B is a block diagram showing one embodiment of portions of a voice operated exchange (VOX) circuit included in the single ear-level device of FIG. 1.</p> <p>Victorian '313 at [0044] Ear-level device 110 includes an external microphone 420, an ear canal microphone 440, an ear canal</p>

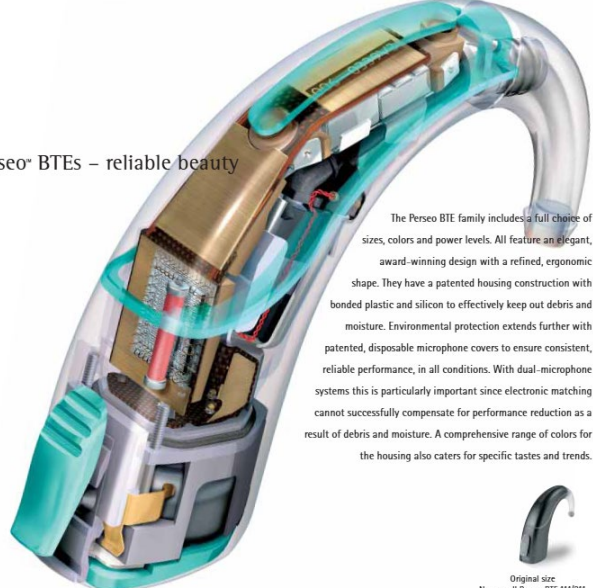
Claims	Prior Art
	<p>speaker 470, a processor 460, a wireless transceiver 450, an antenna 455, a voice operated exchange (VOX) circuit 430, and a battery 480. External microphone 420 picks up ambient sound such that a hearing impaired person wearing ear-level device 110 is not “isolated” when communicating to remote device 190. Thus, ear-level device enables the hearing impaired person to talk to another person either directly (if the other person is nearby) or through a device such as a cell phone. Ear canal microphone 440 detects sound from the ear canal. The sound includes primarily speech of the person wearing ear-level device 110. Ear canal speaker 470 transmits sound received from remote device 190 and/or external microphone 420 to the same ear canal. In one embodiment, ear canal microphone 440 and ear canal speaker 470 are implemented as one physical device. Processor 460 converts the sound picked up from the ear canal to an electrical signal to be transmitted to remote device 190, and converts the signals received from remote device 190 and/or external microphone 420 to a sound audible to the person wearing ear-level device 110. In one embodiment, processor 460 includes one or more of amplification circuitry, filtering circuitry, acoustic feedback reduction circuitry, noise reduction circuitry, and tone control circuitry, among other circuits performing signal processing functions as known in the art. In one embodiment, processor 460 includes a speech recognition module to enhance the audio signal received by and/or transmitted from ear-level device 110. Wireless transceiver 450 and antenna 455 form a telemetry interface supporting telemetry link 195 between ear-level device 110 and remote device 190. Wireless transceiver 450 includes a wireless transmitter and a wireless receiver. The wireless transmitter receives the electrical signal representing the sound picked up from the ear canal from processor 460 and transmits an outgoing wireless signal representing the same sound to remote device 190. The wireless receiver receives an incoming wireless signal representing a remote sound from remote device 190 and transmits an electrical signal representing the remote sound to processor 460. Telemetry link 195 provides for bi-directional communication allowing simultaneous signal transmission, in both directions, between ear-level device 110 and remote device 190. VOX circuit 430 activates a major portion of ear-level device 110 only when sound is detected by at least one of external microphone 420, ear canal microphone 440, and wireless transceiver 450. In one embodiment, VOX circuit 430 includes a voice-controlled switch connecting between battery 480 and the portions of ear-level device 110 that are voice activated. The switch is driven by the sound detected by at least one of external microphone 420, ear canal microphone 440, and wireless transceiver 450. It is turned on when sound detected by external microphone 420 exceeds a predetermined threshold, when sound detected by ear canal microphone 440 exceeds another predetermined threshold, and/or when wireless transceiver 450 receives the incoming wireless signal representing the remote sound. In one embodiment, as discussed with respect to FIG. 4B below, VOX circuit 430 gates or attenuates one or more sounds detected by external microphone 420, ear canal microphone 440, and wireless transceiver 450, to eliminate or reduce echo and ringing caused by the loop formed by ear canal speaker 470 and ear canal microphone 440. Battery 480 supplies the power needed</p>

Claims	Prior Art
	<p>for the operation of ear-level device 110. In one embodiment, battery 480 is a rechargeable battery that can be recharged without being taken out from ear-level device 110.</p> <p>Victorian '313 at [0045]</p> <p>FIG. 4B is a block diagram showing one embodiment of portions of VOX circuit 430. In this embodiment, VOX circuit 430 includes an ambient sound gating module 421 that blanks or attenuates the sound detected by external microphone 420, an occluded sound gating module 441 that blanks or attenuates the sound detected by ear canal microphone 440, a remote sound gating module 451 that blanks or attenuates the sound detected by wireless transceiver 450, and controller 431 to coordinate the states of the gating modules. Because ear canal microphone 440 is capable of detecting the sound delivered to the ear canal from ear canal speaker 470, there is a need to avoid echo and/or ringing. For example, a remote sound from remote device 190 delivered to the ear canal may be detected as an occluded sound and transmitted back to remote device 190, causing echo. The gating modules function in a coordinated manner to eliminate or reduce such echo and/or ringing caused by breaking, or applying resistance into, the audio loop formed by ear canal speaker 470 and ear canal microphone 440.</p> <p>Victorian '313 at [0047]</p> <p>Ambient sound gating module 421 includes a voice-driven switch or attenuator driven by the output of a comparator that compares the amplitude of a detected ambient sound with a predetermined threshold for the ambient sound. When the amplitude of the detected ambient sound exceeds the threshold for the ambient sound, and when controller 431 permits, ambient sound gating module 421 is turned on. When the amplitude of the detected ambient sound drops below the threshold for the ambient sound, ambient sound gating module 421 is turned off. In one embodiment, after the amplitude of the detected ambient sound drops below the threshold for the ambient sound, ambient sound gating module 421 is turned off after a predetermined delay.</p> <p>A POSITA would have been motivated to combine the teachings of the DHA system with Jannard, and Victorian to meet the claim limitations with a reasonable expectation of success for the reasons described above with respect to claims [1a], 26, and 71.</p> <p>To the extent that the DHA system, alone or in view of Jannard, and Victorian, does not disclose any elements of this claim limitation, it would have been obvious to a person having ordinary skill in the art to combine, with the teachings of DHA systems, those elements, as known to them or as disclosed by other prior art identified in the cover pleading.</p>

Claims	Prior Art
<p>[51d] amplify audio output from the first speaker dependent at least in part on the ambient sound level.</p>	<p>DHA systems, alone or in view of Victorian, expressly or inherently discloses “amplify audio output from the first speaker dependent at least in part on the ambient sound level.”</p> <p>Phonak – Perseo Brochure at Page 7</p> <div data-bbox="636 321 737 418" style="display: inline-block; width: 40px; height: 40px; background-color: #FFD700; margin-right: 10px;"></div> <div data-bbox="758 321 1207 386" style="display: inline-block; vertical-align: top;"> <p>Convenience and Comfort Perseo is easier to wear and easier to fit</p> </div> <div data-bbox="1230 321 1331 418" style="display: inline-block; width: 40px; height: 40px; background-color: #0070C0; margin-right: 10px;"></div> <div data-bbox="1352 321 1816 386" style="display: inline-block; vertical-align: top;"> <p>Clarity Enhanced performance in multiple listening situations</p> </div> <div data-bbox="737 443 1816 841" style="display: flex; flex-wrap: wrap;"> <div style="width: 50%;"> <ul style="list-style-type: none"> ■ For true “wear and hear” convenience and comfort in all situations, Perseo automatic operation is managed by Phonak-exclusive PersonalLogic. ■ This breakthrough intelligence for automatic hearing instruments responds to acoustics, the listening situation and personal listening preferences. ■ Perseo is easily programmed for individual needs using the Listening Situation Manager. </div> <div style="width: 50%;"> <ul style="list-style-type: none"> ■ The Clarity Component Enhancement system realized in Digital Perception Processing² (DPP²) increases the clarity of speech and environmental sounds. ■ Rapid-tracking Adaptive digital AudioZoom (dAZ) is automatically activated to improve speech understanding in noisy situations. ■ For extra comfort and ease of listening, the Fine-scale Noise Canceler selectively reduces background noise in 20 channels. ■ Optional MicroLink FM provides the ultimate Speech-to-Noise Ratio enhancement. </div> </div>

Claims	Prior Art
	<p data-bbox="583 175 1062 207">Phonak – Perseo Brochure at Page 13</p> <p data-bbox="638 227 1167 266">Digital Perception Processing² (DPP²)</p>  <p data-bbox="638 646 756 662">Personal Audiogram</p> <p data-bbox="638 727 798 795">LPP Loudness Perception Profile</p> <p data-bbox="865 328 1033 360">Cochlear Perception Model Calculate normal perception</p> <p data-bbox="865 376 1050 560"></p> <p data-bbox="865 587 1050 771"></p> <p data-bbox="865 787 1041 820">Personal Data Calculate personal perception</p> <p data-bbox="1117 328 1327 409">Clarity Component Enhancement: - Identify and enhance Clarity components - Optimize personal perception</p> <p data-bbox="1117 544 1306 820">  </p> <p data-bbox="1369 544 1453 576">Result</p> <p data-bbox="1369 587 1558 771"></p> <p data-bbox="1369 787 1495 820">Optimal cochlea stimulation</p>

Claims	Prior Art
	<p data-bbox="583 175 1064 207">Phonak – Perseo Brochure at Page 14</p> <p data-bbox="632 228 1205 264">Enhanced comfort in all situations</p> <p data-bbox="722 274 1073 297">Fine-scale Noise Canceler in 20 channels</p> <div data-bbox="722 321 1073 459"> <p>Finely scaled frequency resolution in 20 channels gives the Fine-scale Noise Canceler (FNC) the most precise background noise reduction capability. Noise reduction depends on the relative amounts of speech and noise in each channel (graph 1) and the importance of each</p> </div> <div data-bbox="1098 321 1449 459"> <p>channel for speech understanding. The FNC in Perseo ensures maximum comfort and ease of listening. Annoying background sounds are reduced for comfort while the important sounds of speech are preserved for stress-free listening. The maximum amount of noise</p> </div> <div data-bbox="1474 321 1824 459"> <p>reduction can be fitter-programmed to match the individual wearer's tastes and requirements. Thanks to the dual-base automatic, separate settings can be made for relatively quiet and noisy situations, ensuring optimal comfort in multiple situations (graph 2).</p> </div> <div data-bbox="814 508 1824 813"> <p>The figure contains two graphs, Graph 1 and Graph 2, illustrating noise reduction in 20 channels.</p> <p>Graph 1: This graph shows the input signals and the resulting noise reduction. On the left, under the heading "Input", there are two waveforms: a red one labeled "Noise" and a blue one labeled "Speech". To the right of these is a combined waveform. Below this, under the heading "Noise reduction", is a bar chart where the height of each bar represents the level of noise reduction for a specific frequency. The x-axis is labeled "Frequency".</p> <p>Graph 2: This graph shows the noise reduction levels across frequency. It features a shaded area that represents the range of noise reduction. The area is divided into three horizontal sections labeled "Light", "Medium", and "Strong" on the right side. The x-axis is labeled "Frequency". Three red arrows point downwards from the top of the shaded area, indicating a reduction in noise level across the frequency spectrum.</p> </div>

Claims	Prior Art
	<p data-bbox="583 175 1062 207">Phonak – Perseo Brochure at Page 23</p> <div data-bbox="625 224 1260 824">  <p data-bbox="636 383 905 407">Perseo® BTEs – reliable beauty</p> <p data-bbox="968 443 1255 695">The Perseo BTE family includes a full choice of sizes, colors and power levels. All feature an elegant, award-winning design with a refined, ergonomic shape. They have a patented housing construction with bonded plastic and silicon to effectively keep out debris and moisture. Environmental protection extends further with patented, disposable microphone covers to ensure consistent, reliable performance, in all conditions. With dual-microphone systems this is particularly important since electronic matching cannot successfully compensate for performance reduction as a result of debris and moisture. A comprehensive range of colors for the housing also caters for specific tastes and trends.</p> <p data-bbox="1140 800 1247 824">Original size New small Perseo BTE 111/211</p> </div> <p data-bbox="583 846 1625 878">Valente - Effect of type of noise and loudspeaker array on the performance at 399</p> <p data-bbox="636 898 1022 1157">Recently, Phonak introduced its digital hearing aid in which one feature is an adaptive directional microphone. This adaptive microphone reportedly ensures identical sensitivity and phase characteristics of the dual microphones and a mechanism to <i>automatically</i> change the polar pattern from omnidirectional to any directional pattern while compensating for the changes in low-frequency response associated with each pattern. This hearing aid can be programmed</p>

Claims

Prior Art

Phonak – The Hearing Review Article re Perseo (Sep. 5, 2003) at Page 2

Perseo is a 20-channel digital hearing instrument. The default base-program is optimized for speech understanding in quiet environments. A second base-program with different frequency response settings and optional adaptive directional microphones provides additional support in noisy situations. The acoustic environment is continuously monitored and the appropriate base-program is automatically selected. This dual-base automatic function allows for independent fine tuning in one listening situation without compromising performance in the other. The automatic performance is designed to be easily personalized for the individual needs and preferences of the user.

Phonak – Collections 2002 Product List at Page 6

Elite



Digital Perception Processing (DPP)



AutoSelect (AS)

<p>New Phonak sound: unmatched clarity Digital Perception Processing (DPP) is based on a cochlear sound analysis across 20 critical bands - exactly like the analysis performed by the healthy human ear. But Clari gives a step further than previous overlapping multi-band filter systems. DPP models the active, load-free dependent cochlear signal processing in all 20 critical bands.</p>	<p>20 critical bands: overlapping and dynamic The cochlear filter is overlapping, but the degree of overlap is dependent on the frequency, level, and bandwidth of the sound signal. This is precisely how Clari's filter work. Digital Perception Processing computes the cochlear excitation patterns and thus determines sound perception. DPP is the first to mimic the healthy human ear in this way.</p>	<p>Automatic, intelligent program selection Clari listens and automatically activates the optimal DPP algorithm for the acoustic environment: Program 1 with QuietAdapt for quiet situations or Program 2 with NoiseAdapt and noise control for noisy environments.</p> <p>In quiet situations QuietAdapt DPP is designed for clear hearing in quiet surroundings. QuietAdapt controls gain in 20 critical bands to correct for loudness and provide a natural sound quality. By preserving the perception patterns of complex sounds, QuietAdapt assures that low-level signal components are not over-</p>	<p>amplified to the point where they would adversely affect speech intelligibility. This allows Clari to deliver only those sound components that would be naturally audible.</p> <p>In noise NoiseAdapt DPP is the ideal program for maximizing audibility in noisy environments. The special NoiseAdapt algorithm works with the fine-scale Noise Canceller - and if available with the Adaptive digital AudioZoom feature - to optimize speech understanding in noise.</p>
--	---	--	--



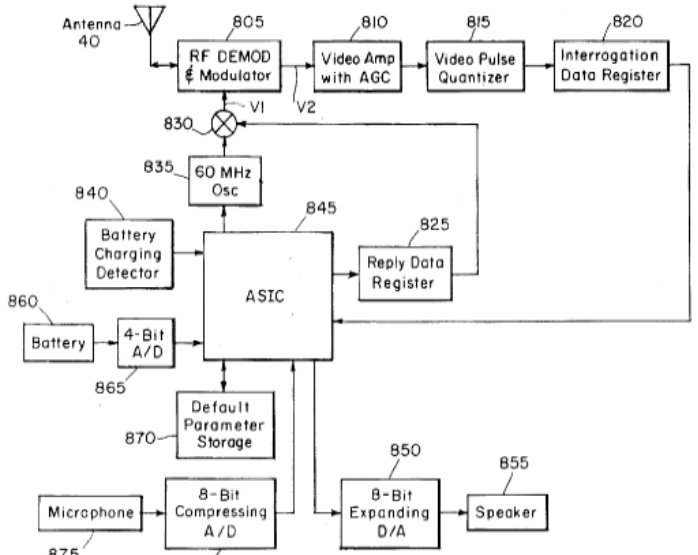
Oticon – DigiFocus Description at Page 1

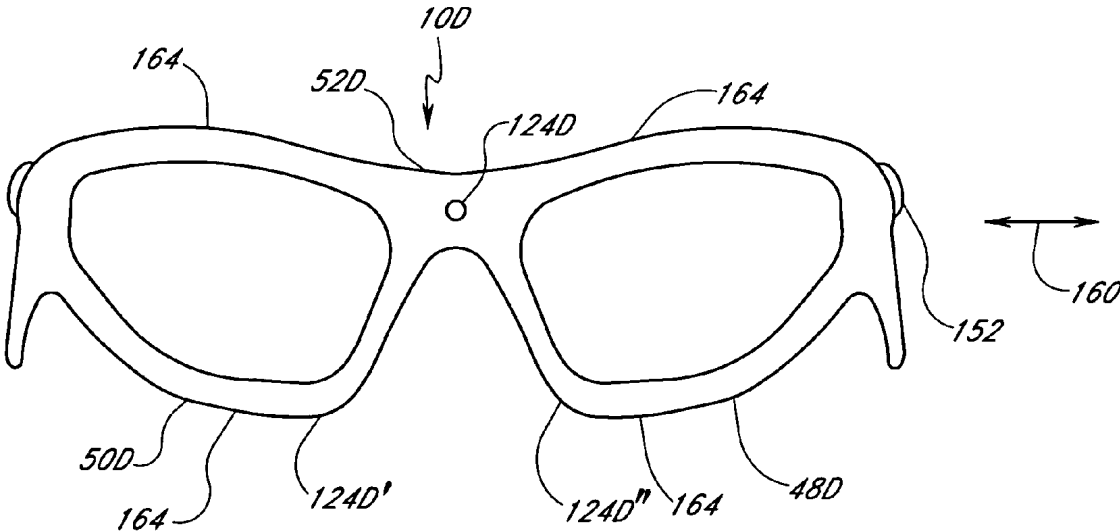
Inside each DigiFocus, the digital chips have the same processing power as a fast laptop computer. This power is used to analyse the incoming sound signal, and speech is extracted through a unique system called ASA - Adaptive Speech Alignment. ASA always makes speech audible, amplifying it as much as necessary without making it uncomfortably loud or distorted. Sounds are reproduced very naturally - making speech far easier to understand than with any other hearing instrument on the market.

Claims	Prior Art
	<p data-bbox="583 175 1976 321">Oticon – LA Time Article re THE GOODS: The Sound Wave of the Future (Nov. 10 1995) at 1 “Today’s amplifiers are not smart,” said Peter Hahn, managing director of Oticon’s U.S. headquarters in Somerset, N.J. “The digital processor can organize frequencies and manipulate speech signals to separate them from background noise.”</p> <p data-bbox="583 375 1976 1474">Victorian ’313 at [0044] Ear-level device 110 includes an external microphone 420, an ear canal microphone 440, an ear canal speaker 470, a processor 460, a wireless transceiver 450, an antenna 455, a voice operated exchange (VOX) circuit 430, and a battery 480. External microphone 420 picks up ambient sound such that a hearing impaired person wearing ear-level device 110 is not “isolated” when communicating to remote device 190. Thus, ear-level device enables the hearing impaired person to talk to another person either directly (if the other person is nearby) or through a device such as a cell phone. Ear canal microphone 440 detects sound from the ear canal. The sound includes primarily speech of the person wearing ear-level device 110. Ear canal speaker 470 transmits sound received from remote device 190 and/or external microphone 420 to the same ear canal. In one embodiment, ear canal microphone 440 and ear canal speaker 470 are implemented as one physical device. Processor 460 converts the sound picked up from the ear canal to an electrical signal to be transmitted to remote device 190, and converts the signals received from remote device 190 and/or external microphone 420 to a sound audible to the person wearing ear-level device 110. In one embodiment, processor 460 includes one or more of amplification circuitry, filtering circuitry, acoustic feedback reduction circuitry, noise reduction circuitry, and tone control circuitry, among other circuits performing signal processing functions as known in the art. In one embodiment, processor 460 includes a speech recognition module to enhance the audio signal received by and/or transmitted from ear-level device 110. Wireless transceiver 450 and antenna 455 form a telemetry interface supporting telemetry link 195 between ear-level device 110 and remote device 190. Wireless transceiver 450 includes a wireless transmitter and a wireless receiver. The wireless transmitter receives the electrical signal representing the sound picked up from the ear canal from processor 460 and transmits an outgoing wireless signal representing the same sound to remote device 190. The wireless receiver receives an incoming wireless signal representing a remote sound from remote device 190 and transmits an electrical signal representing the remote sound to processor 460. Telemetry link 195 provides for bi-directional communication allowing simultaneous signal transmission, in both directions, between ear-level device 110 and remote device 190. VOX circuit 430 activates a major portion of ear-level device 110 only when sound is detected by at least one of external microphone 420, ear canal microphone 440, and wireless transceiver 450. In one embodiment, VOX circuit 430 includes a voice-controlled switch connecting between battery 480 and the portions of ear-level device 110 that are voice activated. The switch is driven by the sound detected by at least one of external microphone 420, ear canal microphone</p>

Claims	Prior Art
	<p>440, and wireless transceiver 450. It is turned on when sound detected by external microphone 420 exceeds a predetermined threshold, when sound detected by ear canal microphone 440 exceeds another predetermined threshold, and/or when wireless transceiver 450 receives the incoming wireless signal representing the remote sound. In one embodiment, as discussed with respect to FIG. 4B below, VOX circuit 430 gates or attenuates one or more sounds detected by external microphone 420, ear canal microphone 440, and wireless transceiver 450, to eliminate or reduce echo and ringing caused by the loop formed by ear canal speaker 470 and ear canal microphone 440. Battery 480 supplies the power needed for the operation of ear-level device 110. In one embodiment, battery 480 is a rechargeable battery that can be recharged without being taken out from ear-level device 110.</p> <p>Victorian '313 at [0045]</p> <p>FIG. 4B is a block diagram showing one embodiment of portions of VOX circuit 430. In this embodiment, VOX circuit 430 includes an ambient sound gating module 421 that blanks or attenuates the sound detected by external microphone 420, an occluded sound gating module 441 that blanks or attenuates the sound detected by ear canal microphone 440, a remote sound gating module 451 that blanks or attenuates the sound detected by wireless transceiver 450, and controller 431 to coordinate the states of the gating modules. Because ear canal microphone 440 is capable of detecting the sound delivered to the ear canal from ear canal speaker 470, there is a need to avoid echo and/or ringing. For example, a remote sound from remote device 190 delivered to the ear canal may be detected as an occluded sound and transmitted back to remote device 190, causing echo. The gating modules function in a coordinated manner to eliminate or reduce such echo and/or ringing caused by breaking, or applying resistance into, the audio loop formed by ear canal speaker 470 and ear canal microphone 440.</p> <p>Victorian '313 at [0047]</p> <p>Ambient sound gating module 421 includes a voice-driven switch or attenuator driven by the output of a comparator that compares the amplitude of a detected ambient sound with a predetermined threshold for the ambient sound. When the amplitude of the detected ambient sound exceeds the threshold for the ambient sound, and when controller 431 permits, ambient sound gating module 421 is turned on. When the amplitude of the detected ambient sound drops below the threshold for the ambient sound, ambient sound gating module 421 is turned off. In one embodiment, after the amplitude of the detected ambient sound drops below the threshold for the ambient sound, ambient sound gating module 421 is turned off after a predetermined delay.</p>

Claims	Prior Art
	<p>A POSITA would have been motivated to combine the teachings of the DHA system with Victorian to meet the claim limitations with a reasonable expectation of success for the reasons described above with respect to claims [1a] and 71.</p> <p>To the extent that the DHA system, alone or in view of Victorian, does not disclose any elements of this claim limitation, it would have been obvious to a person having ordinary skill in the art to combine, with the teachings of DHA systems, those elements, as known to them or as disclosed by other prior art identified in the cover pleading.</p>
<p>[52a] A personal audio system as recited in claim 42, wherein the first head-worn structure of the wearable electronic device comprises: a battery internal to the first head-worn structure; and</p>	<p>DHA systems, alone or in view of Anderson and Jannard, expressly or inherently discloses “[a] personal audio system as recited in claim 42, wherein the first head-worn structure of the wearable electronic device comprises: a battery internal to the first head-worn structure.”</p> <p><i>See above at claim [25a].</i></p> <p>A POSITA would have been motivated to combine the teachings of the DHA system with Anderson and Jannard to meet the claim limitations with a reasonable expectation of success for the reasons described above with respect to claim [1a] and [26a].</p> <p>To the extent that the DHA system, alone or in view of Anderson and Abreu, does not disclose any elements of this claim limitation, it would have been obvious to a person having ordinary skill in the art to combine, with the teachings of DHA systems, those elements, as known to them or as disclosed by other prior art identified in the cover pleading.</p>
<p>[52b] at least one conductive element that facilitates electrical coupling with the battery, the at least one conductive element is at least partially internal to the first head-worn structure but remains at least partially exposed external to the first head-worn structure, and the at least one conductive element comprises at least one conductive pad,</p>	<p>DHA systems, alone or in view of Anderson and Jannard, expressly or inherently discloses “at least one conductive element that facilitates electrical coupling with the battery, the at least one conductive element is at least partially internal to the first head-worn structure but remains at least partially exposed external to the first head-worn structure, and the at least one conductive element comprises at least one conductive pad.”</p>

Claims	Prior Art
	<p data-bbox="583 175 898 207">Anderson '783 at FIG. 8</p>  <p data-bbox="945 795 1050 828">FIG. 8</p> <p data-bbox="583 844 919 876">Anderson '783 at 3:40–41</p> <p data-bbox="630 885 1963 950">FIG. 8 is a block diagram showing details of an earpiece for a full-featured preferred embodiment of the invention.</p> <p data-bbox="583 1006 919 1039">Anderson '783 at 6:55–57</p> <p data-bbox="630 1047 1932 1112">Earpiece power consumption is preferably minimized to extend battery life, and use of a rechargeable battery is desirable.</p>

Claims	Prior Art
	<p data-bbox="583 175 884 207">Jannard '526 at FIG. 11</p>  <p data-bbox="978 769 1157 826">FIG. 11</p> <p data-bbox="583 846 894 878">Jannard '526 at 4:27–28</p> <p data-bbox="632 883 1604 915">FIG. 11 is a front elevational view of the audio device illustrated in FIG. 10.</p> <p data-bbox="583 971 936 1003">Jannard '526 at 24:56–25:3</p> <p data-bbox="632 1008 1976 1300">The audio device 10D can also include power contacts 164 for recharging any rechargeable batteries connected thereto. For example, the power contacts 164 can be disposed on a lower edge of the orbitals 48D, 50D. Thus, with an appropriate recharging cradle (not shown), the audio device 10D can be laid on the cradle, thereby making contact between the power contacts 164 and corresponding contacts in the cradle (not shown). Alternatively, power contacts can be provided in numerous other locations as desired. For example, the power contacts 164 can be disposed at the ends of the ear stems 54D, 56D. A corresponding cradle can include two vertically oriented holes into which the ear stems are inserted for recharging. In this configuration, the lens within the orbitals 48D, 50D would face directly upwardly.</p> <p data-bbox="583 1354 1986 1459">A POSITA would have been motivated to combine the teachings of the DHA system with Anderson and Jannard to meet the claim limitations with a reasonable expectation of success for the reasons described above with respect to claim [1a] and [26a].</p>

Claims	Prior Art
	<p>To the extent that the DHA system, alone or in view of Anderson and Abreu, does not disclose any elements of this claim limitation, it would have been obvious to a person having ordinary skill in the art to combine, with the teachings of DHA systems, those elements, as known to them or as disclosed by other prior art identified in the cover pleading.</p>
<p>[52c] wherein the battery is able to be electrically coupled via the at least one conductive element to a charging apparatus and charged when at least the first head-worn structure of the head-worn wearable electronic device is placed within the charging apparatus, and</p>	<p>DHA systems, alone or in view of Anderson and Jannard, expressly or inherently discloses “wherein the battery is able to be electrically coupled via the at least one conductive element to a charging apparatus and charged when at least the first head-worn structure of the head-worn wearable electronic device is placed within the charging apparatus.”</p> <p><i>See above at claim [26b].</i></p> <p>A POSITA would have been motivated to combine the teachings of the DHA system with Anderson and Jannard to meet the claim limitations with a reasonable expectation of success for the reasons described above with respect to claim [1a] and [26a].</p> <p>To the extent that the DHA system, alone or in view of Anderson and Abreu, does not disclose any elements of this claim limitation, it would have been obvious to a person having ordinary skill in the art to combine, with the teachings of DHA systems, those elements, as known to them or as disclosed by other prior art identified in the cover pleading.</p>
<p>[52d] wherein the electrical coupling of the at least one conductive element with the charging apparatus is supported by at least a magnetic element.</p>	<p>DHA systems, alone or in view of Anderson and Jannard, expressly or inherently discloses “wherein the electrical coupling of the at least one conductive element with the charging apparatus is supported by at least a magnetic element.”</p> <p><i>See above at claim [26c].</i></p> <p>A POSITA would have been motivated to combine the teachings of the DHA system with Anderson and Jannard to meet the claim limitations with a reasonable expectation of success for the reasons described above with respect to claim [1a] and [26a].</p> <p>To the extent that the DHA system, alone or in view of Anderson and Jannard , does not disclose any elements of this claim limitation, it would have been obvious to a person having ordinary skill in the art to combine, with the teachings of DHA systems, those elements, as known to them or as disclosed by other prior art identified in the cover pleading.</p>

Claims	Prior Art
<p>[53a] A personal audio system as recited in claim 52, wherein the head-worn wearable electronic device comprises: a microphone at least partially internal to the first head-worn structure,</p>	<p>DHA systems, alone or in view of Anderson, expressly or inherently discloses “[a] personal audio system as recited in claim 52, wherein the head-worn wearable electronic device comprises: a microphone at least partially internal to the first head-worn structure.”</p> <p><i>See above at claim 6.</i></p> <p>A POSITA would have been motivated to combine the teachings of the DHA system with Anderson to meet the claim limitations with a reasonable expectation of success for the reasons described above with respect to claim [1a] and [26a].</p> <p>To the extent that the DHA system, alone or in view of Anderson, does not disclose any elements of this claim limitation, it would have been obvious to a person having ordinary skill in the art to combine, with the teachings of DHA systems, those elements, as known to them or as disclosed by other prior art identified in the cover pleading.</p>
<p>[53b] wherein the head-worn wearable electronic device is configured to receive a voice input via the microphone, determine a requested operation from the voice input using voice recognition, and control an operation by the head-worn wearable electronic device based on the requested operation.</p>	<p>DHA systems, alone or in view of Anderson and Jannard or Seshadri, expressly or inherently discloses “wherein the head-worn wearable electronic device is configured to receive a voice input via the microphone, determine a requested operation from the voice input using voice recognition, and control an operation by the head-worn wearable electronic device based on the requested operation.”</p> <p>Anderson ’783 at 22:11–61</p> <p>Stand-alone earpiece operation: If an earpiece cannot communicate with an RPU for any reason, stand-alone earpiece operation as a CIC hearing aid is generally desirable. Examples include occasions where the RPU must be turned off by pushing the RPU power-off key 78 (e.g., situations where operation of wireless transmitters is prohibited such as aircraft takeoffs and classified military briefings), RPU battery or other RPU electronics failure, misplacement of an RPU, and primary wireless link 17 interference. Note that a microphone located in the earpiece 10 is preferred to a body-worn microphone, since only an earpiece microphone 12 can be used to provide the desired default CIC hearing aid function (earpiece microphones also provide superior audio localization capabilities as the microphone moves with the user's head, and provide reduced wind and clothing noise compared with body-worn microphones). Parameters that control the stand-alone earpiece operation are stored in the earpiece CIC default parameter storage 870, which is accessed by the ASIC 845 to implement the desired stand-alone functions as described earlier. CIC default parameters include, but are not limited to, settings for the nominal and maximum allowable earpiece volume. Other information contained in the default parameter storage 870 includes, but is not limited to, the communication address and L/R (left/right) selection bit. To ensure that the CIC default parameter storage information cannot be changed inadvertently, changes are only allowed when the earpiece is not in use as, for example, when the earpiece battery 860 is being charged. A battery charging detector 840, that is sensitive to the high-level magnetic induction fields commonly used for</p>

Claims	Prior Art
	<p>earpiece battery charging, indicates to the ASIC 845 that the default parameter storage 870 values may be changed upon receipt of an appropriate command from the RPU. Specifically, to change information in the earpiece CIC default parameter storage 870, a user first activates the earpiece battery charging system, e.g., by placing the earpiece in a recharger/storage tray well known in the art. The RPU keyboard 946 is then used to enter a command to the DSP 948 that indicates the user's desire to change information in the default parameter storage 870, followed by the new information to be stored. The DSP 948 causes the RPU to transmit an interrogation with a reserved address code to which all earpieces being charged within the communication range of the RPU are receptive. The information that is to ultimately be loaded into the earpiece default parameter storage 870 is then transmitted by the RPU using the auxiliary bit field of many successive interrogations, said interrogations are received by the earpiece, and the desired information is loaded into the CIC default parameter storage 870 under control of the ASIC 845.</p> <p>Anderson '783 at 25:15–27</p> <p>To obtain solicited information, user requests are processed in a variety of means well known in the art, including inputs via one or more pushbuttons on the RPU and voice recognition computer programs implemented in the RPU electronics. Voice recognition can be used almost exclusively, if desired, thereby eliminating nearly all RPU pushbutton controls. Voice commands recognized by the RPU can cause the RPU to initiate actions such as changing hearing aid parameter settings (e.g., volume), or activate appliances and environmental controls (e.g., aids for the physically handicapped) in a distributed wireless network via the secondary wireless link.</p> <p>Anderson '783 at 25:29–30</p> <p>The verbal commands are received as audio signals by the earpiece 10 and transmitted to the RPU 16.</p> <p>Jannard '526 at 25:53–26:20</p> <p>In another embodiment, a further modification of the audio devices 10, 10A, 10B, 10C, and 10D is provided wherein the audio devices include at least two banks of microphones, with one bank acting as a speaker of received and one bank providing an ambient noise-cancellation function. The microphone banks can be positioned at any suitable location or combination of locations (e.g., on the audio device, within the audio device, opposing sides of the audio device, or the like). In one embodiment, automatic switching of the speaking-microphone and noise-canceling-microphone banks' functions advantageously enhances ease of use. In a further embodiment, the microphone banks can be arranged in an array to be used in conjunction with algorithms to discern, reduce, and/or eliminate noise for the purpose of voice recognition. For example, in one embodiment, such microphone banks can include ASIC-based noise-</p>

Claims	Prior Art
	<p>canceling technology, such as is available in chips from Andrea Electronics Corporation (AEC), to enable voice recognition in ambient noise up to about 130 Db or more. In another embodiment, microphone banks can be arranged in any suitable combination of linear or non-linear arrays to be used in conjunction with algorithms to discern, reduce, and/or eliminate noise in any suitable manner. In another embodiment, audio/proximity sensors can advantageously trigger the appropriate functionality in a specific bank. In another embodiment, a noise-canceling microphone can be provided in connection with a cord or other microphones described above. For example, without limitation, a series of miniature microphones can be supported down a cord from the audio device, separated by desired distances, and aimed in different directions. In another implementation, one or more of the microphones can be for verbal input from the user, and one or more others of the microphones, or the same microphone, can also be for noise-cancellation purposes.</p> <p>Jannard '526 at 30:53–31:2</p> <p>With reference to FIG. 9A, the interface electronics 116 can include audio or “aural” menus that can be selected by user. For example, a user can initiate an audio menu by depressing the button 150 (FIGS. 10–12). Upon initiation of the audio menus, the interface electronics 116 can send an audio signal to the speakers 14, 14A, 14B, 14C, 14D, 14E, 16, 16A, 16B, 16C, 16D, 16E including a humanoid voice. The voice signal can indicate that a first menu option is available. For example, but without limitation, the first menu choice can be to initiate a call. Thus, when the user pushes the button 150 the first time, the user will hear the words “initiate a call,” emanating from the speakers 14, 14A, 14B, 14C, 14D, 14E, 16, 16A, 16B, 16C, 16D, 16E. If the user wishes to initiate a call, the user can depress the button 150 again which will send the appropriate AT command to the transceiver 114 so as to transmit the proper AT code to the cellular phone source device S, B (FIG. 8).</p> <p>Jannard '526 at 32:4–13</p> <p>For example, if a user chooses a voice command option, the user could issue commands such as, for example, but without limitation, “phonebook” or “call alpha.” With a source device such as a cellular phone, that has a speech recognition engine and that is properly trained to recognize the voice of the user, the user can automatically enter the phonebook mode or directly call the phonebook listing “alpha,” of course, as is apparent to one of ordinary skill in the art, such a voice command protocol could be used to issue other commands as well.</p> <p>Jannard '526 at 32:14–19</p> <p>In another alternative, the interface electronics 116 can include a speech recognition engine and audio</p>

Claims	Prior Art
	<p>menus. In this alternative, the interface electronics 116 can recognize speech from the user, convert the speech to AT commands, and control this source device using a standard AT command protocol.</p> <p>Seshadri '839 at [0042] Input module 128, which may be a keypad, touch screen, voice recognition circuit, or other like user interfaces, receives user commands and produces digital command messages 124 there from. Such digital command messages 124 includes, but are not limited to, packet size, synchronization information, frequency hopping initiation information, timeslot allocation information, link establishment information, piconet address information, fast-forward, play, pause, volume adjust, record, stop and rewind.</p> <p>Seshadri '839 at [0047] Detachable ear-piece 156 and microphone 154 may have on-chip operations to support call conferencing, call waiting, flash, and other features associated with telephones. These functions may be accessed and reviewed by a user interface 158 and display 160 within host device 152 or a user interface and display located on either detachable ear-piece 156 or microphone 154. The user interface and display, located on either the host device or detachable ear-piece 156 and microphone 154 may have a display and button(s) that may be used to program device, perform directory functions including selecting number to call, view caller ID, initiate call waiting, or initiate call conferencing. Additionally, circuitry within the ear-piece 156 and microphone 154 may enable voice activated dialing. The actual voice recognition could be performed within ear-piece 156, microphone 154, or host device 152. Thus, ear-piece 156 and microphone 154 may act to initiate calls and receive calls.</p> <p>A POSITA would have been motivated to combine the teachings of the DHA system with Anderson and Jannard or Seshadri to meet the claim limitations with a reasonable expectation of success for the reasons described above with respect to claim [1a] and [26a].</p> <p>To the extent that the DHA system, alone or in view of Anderson and Jannard or Seshadri, does not disclose any elements of this claim limitation, it would have been obvious to a person having ordinary skill in the art to combine, with the teachings of DHA systems, those elements, as known to them or as disclosed by other prior art identified in the cover pleading.</p>
[54a] A personal audio system as recited in claim 42, wherein the first head-worn structure comprises: an infrared emitter	DHA systems, alone or in view of Anderson and Abreu or Irvin or Luo, expressly or inherently discloses “[a] personal audio system as recited in claim 42, wherein the first head-worn structure comprises: an infrared

Claims	Prior Art
<p>configured to at least direct infrared signals towards the user when at least the first head-worn structure of the head-worn wearable electronic device is being worn; and</p>	<p>emitter configured to at least direct infrared signals towards the user when at least the first head-worn structure of the head-worn wearable electronic device is being worn.”</p> <p>Abreu '960 at Abstract</p> <p>The support structure includes patches, clips, eyeglasses, head mounted gear and the like, containing passive or active sensors positioned at the end of the tunnel with sensing systems positioned on and accessing a physiologic tunnel.</p> <p>Abreu '960 at 2:62–3:12</p> <p>Included in the present invention are support structures for positioning a sensor on the skin at the end of the tunnel. The present invention discloses devices directed at measuring brain temperature, brain function, metabolic function, hydrodynamic function, hydration status, hemodynamic function, body chemistry and the like. The components include devices and methods for evaluating biological parameters using patches, clips, eyeglasses, head mounted gear and the like with sensing systems adapted to access physiologic tunnels to provide precise and clinically useful information about the physiologic status of the wearer and for enhancing the safety and performance of said wearer, and helping to enhance and preserve the life of said wearer by providing adequate reporting means and alert means relating to the biological parameter being monitored. Other components provide for producing direct or indirect actions, acting on another device, or adjusting another device or article of manufacture based on the biological parameter measured.</p> <p>Abreu '960 at 44:48–51</p> <p>The support structure includes at least one radiation source from infrared to visible light which interacts with the substance being measured at the BTE tunnel and a detector for collecting the resulting radiation.</p> <p>Abreu '960 at 44:61–67</p> <p>The present invention includes means for directing preferably near-infrared energy into the surface of the skin at the end of the BTE tunnel, means for analyzing and converting the reflectance or back scattered spectrun into the concentration of the substance measured and support means for positioning the light source and detector means adjacent to the surface of the skin at the BTE tunnel site.</p> <p>Irvin '332 at 6:3–33</p> <p>FIG. 4C illustrates the power control circuit 210 coupled with a contact switch functioning as the sensor 220. The contact switch includes a first and second contact connected to the power control circuit 210,</p>

Claims	Prior Art
	<p>and connecting member 222 actuated when the user places the wireless headset 100 on their body. When actuated, the connecting member 222 places the first and second contacts of sensor 220 in electrical connection. Ideally, the sensor 220, when configured as a contact sensor, is placed in the earpiece 130, in the retaining member 140, or in the retaining pad 150 shown in FIG. 1B. Preferably, when sensor 220 is configured as a contact-based proximity sensor, the contacting element of sensor 220 is sealed in a manner that minimizes contaminant ingress. In some embodiments, the contact sensor may be actuated by the user resiliently expanding the retaining member 140 or other element of the wireless headset 100 for mounting upon their body.</p> <p>FIG. 4D illustrates the power control circuit 210 coupled with a proximity sensor 220 not based on mechanical actuation. The proximity sensor 220 may be an active or passive IR sensor, may be a resistive or other type of thermal sensor, or may be an acoustic sensor. As a proximity sensor, sensor 220 provides indication that a user has placed the wireless headset 100 upon their body. As such, any of the foregoing or other presence sensing techniques may be used to implement sensor 220. Preferably, when configured as a proximity sensor, sensor 220 requires minimal power and provides a sensor output signal based on proximity to the user's body. Depending upon the sensing technology chosen, the control circuit 210 may be required to periodically stimulate sensor 220 to generate a detection response.</p> <p>Luo '446 at [0036]</p> <p>The invention generally relies on shining light having a particular wavelength on human skin as well as the reflectance properties of human skin. The surface reflection of the human skin occurs at the surface of the epidermis and is approximately independent of the lighting wavelength and independent of human race. In the infrared (IR) range the wavelengths that can be used are greater than 800 nano-meters and the surface reflectance of human skin is more or less constant and close to 50%. The same hold true for orange and red light in the visible light range, namely wavelengths between approximately 600 and 800 nano-meters can be used, although the surface reflectance for darker skin is reduced for shorter wavelengths. Surface reflectance for dark skin is approximately 25% at 700 nano-meters. In the embodiments shown herein, the energy of choice in the sensor unit 18 is preferably infrared (IR) energy. However, long wavelength visible light energy can also be used as discussed below.</p> <p>Luo '446 at [0038]</p> <p>The optical window 56 is placed at a certain location on the shell 58 of the hearing aid 10. The location of the optical window 56 depends on whether the hearing aid is a BTE, ITE, ITC or CIC hearing aid. In addition, the location is chosen to minimize the distance between the optical window 56 and the skin 64. For optimal reflection of IR signals back through the optical window, the hearing aid user's skin (i.e. the reflecting surface) is required to be immediately over the optical window 56; otherwise the reflected IR</p>

Claims	Prior Art
	<p>energy will not be reflected back towards the reception unit 62. The optical window 56 is typically a small window having a diameter of approximately 1 mm for example. The optical window 56 can be made from IR grade glass or other suitable material that allows for the passage of infrared energy (a different material would be used if visible light is used rather than infrared light). The optical window 56 should be kept clean at all times in order to prevent the emitted infrared energy from being reflected back due to dirt and the like that may accumulate on the optical window 56 over time. The location of the optical window 56 will be discussed in further detail below.</p> <p>A POSITA would have been motivated to combine the teachings of the DHA system with Anderson and Abreu or Irvin or Luo to meet the claim limitations with a reasonable expectation of success for the reasons described above with respect to claim [1a] and [23a].</p> <p>To the extent that the DHA system, alone or in view of Anderson and Abreu or Irvin or Luo, does not disclose any elements of this claim limitation, it would have been obvious to a person having ordinary skill in the art to combine, with the teachings of DHA systems, those elements, as known to them or as disclosed by other prior art identified in the cover pleading.</p>
<p>[54b] an infrared detector configured to detect infrared radiation reflected back from the user when at least the first head-worn structure of the head-worn wearable electronic device is being worn.</p>	<p>DHA systems, alone or in view of Anderson and Abreu or Irvin or Luo, expressly or inherently discloses “an infrared detector configured to detect infrared radiation reflected back from the user when at least the first head-worn structure of the head-worn wearable electronic device is being worn.”</p> <p>Abreu ’960 at 44:48–51</p> <p>The support structure includes at least one radiation source from infrared to visible light which interacts with the substance being measured at the BTE tunnel and a detector for collecting the resulting radiation.</p> <p>Abreu ’960 at 44:61–67</p> <p>The present invention includes means for directing preferably near-infrared energy into the surface of the skin at the end of the BTE tunnel, means for analyzing and converting the reflectance or back scattered spectrun into the concentration of the substance measured and support means for positioning the light source and detector means adjacent to the surface of the skin at the BTE tunnel site.</p> <p>Irvin ’332 at 6:3–33</p> <p>FIG. 4C illustrates the power control circuit 210 coupled with a contact switch functioning as the sensor 220. The contact switch includes a first and second contact connected to the power control circuit 210, and connecting member 222 actuated when the user places the wireless headset 100 on their body. When</p>

Claims	Prior Art
	<p>actuated, the connecting member 222 places the first and second contacts of sensor 220 in electrical connection. Ideally, the sensor 220, when configured as a contact sensor, is placed in the earpiece 130, in the retaining member 140, or in the retaining pad 150 shown in FIG. 1B. Preferably, when sensor 220 is configured as a contact-based proximity sensor, the contacting element of sensor 220 is sealed in a manner that minimizes contaminant ingress. In some embodiments, the contact sensor may be actuated by the user resiliently expanding the retaining member 140 or other element of the wireless headset 100 for mounting upon their body.</p> <p>FIG. 4D illustrates the power control circuit 210 coupled with a proximity sensor 220 not based on mechanical actuation. The proximity sensor 220 may be an active or passive IR sensor, may be a resistive or other type of thermal sensor, or may be an acoustic sensor. As a proximity sensor, sensor 220 provides indication that a user has placed the wireless headset 100 upon their body. As such, any of the foregoing or other presence sensing techniques may be used to implement sensor 220. Preferably, when configured as a proximity sensor, sensor 220 requires minimal power and provides a sensor output signal based on proximity to the user's body. Depending upon the sensing technology chosen, the control circuit 210 may be required to periodically stimulate sensor 220 to generate a detection response.</p> <p>Luo '446 at [0036]</p> <p>The invention generally relies on shining light having a particular wavelength on human skin as well as the reflectance properties of human skin. The surface reflection of the human skin occurs at the surface of the epidermis and is approximately independent of the lighting wavelength and independent of human race. In the infrared (IR) range the wavelengths that can be used are greater than 800 nano-meters and the surface reflectance of human skin is more or less constant and close to 50%. The same hold true for orange and red light in the visible light range, namely wavelengths between approximately 600 and 800 nano-meters can be used, although the surface reflectance for darker skin is reduced for shorter wavelengths. Surface reflectance for dark skin is approximately 25% at 700 nano-meters. In the embodiments shown herein, the energy of choice in the sensor unit 18 is preferably infrared (IR) energy. However, long wavelength visible light energy can also be used as discussed below.</p> <p>Luo '446 at [0038]</p> <p>The optical window 56 is placed at a certain location on the shell 58 of the hearing aid 10. The location of the optical window 56 depends on whether the hearing aid is a BTE, ITE, ITC or CIC hearing aid. In addition, the location is chosen to minimize the distance between the optical window 56 and the skin 64. For optimal reflection of IR signals back through the optical window, the hearing aid user's skin (i.e. the reflecting surface) is required to be immediately over the optical window 56; otherwise the reflected IR energy will not be reflected back towards the reception unit 62. The optical window 56 is typically a small</p>

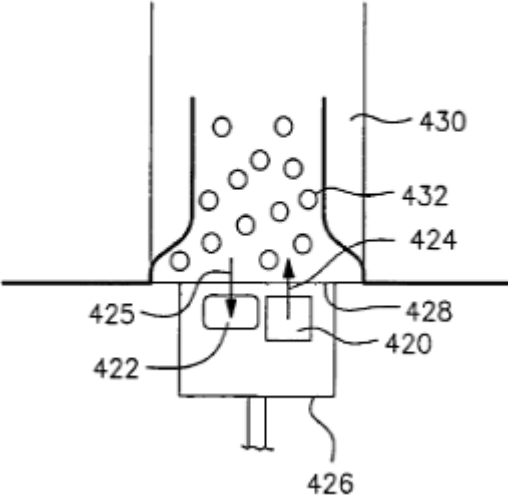
Claims	Prior Art
	<p>window having a diameter of approximately 1 mm for example. The optical window 56 can be made from IR grade glass or other suitable material that allows for the passage of infrared energy (a different material would be used if visible light is used rather than infrared light). The optical window 56 should be kept clean at all times in order to prevent the emitted infrared energy from being reflected back due to dirt and the like that may accumulate on the optical window 56 over time. The location of the optical window 56 will be discussed in further detail below.</p> <p>A POSITA would have been motivated to combine the teachings of the DHA system with Anderson and Abreu or Irvin or Luo to meet the claim limitations with a reasonable expectation of success for the reasons described above with respect to claim [1a] and [23a].</p> <p>To the extent that the DHA system, alone or in view of Anderson and Abreu or Irvin or Luo, does not disclose any elements of this claim limitation, it would have been obvious to a person having ordinary skill in the art to combine, with the teachings of DHA systems, those elements, as known to them or as disclosed by other prior art identified in the cover pleading.</p>
<p>55. A personal audio system as recited in claim 54, wherein electronic circuitry within the first head-worn structure is configured to facilitate determining whether the first head-worn structure is being worn by the user based on at least the detected infrared radiation.</p>	<p>DHA systems, alone or in view of Anderson and Abreu or Irvin or Luo, expressly or inherently discloses “[a] personal audio system as recited in claim 54, wherein electronic circuitry within the first head-worn structure is configured to facilitate determining whether the first head-worn structure is being worn by the user based on at least the detected infrared radiation.”</p> <p>Abreu ’960 at 15:57–60 The eyeglasses can be incorporated with a detecting system to send a signal if the eyeglasses were removed or if the temperature sensor is not capturing signals in a proper manner.</p> <p>Abreu ’960 at 15:60–65 By way of illustration, but not of limitation, pressuse sensing means can be incorporated in the end of the temples to detect if the sunglasses are being worn, and an abrupt drop in the pressure signal indicates glasses were removed or misplacement of the sensor can also generate an identifiable signal.</p> <p>Abreu ’960 at 18:56–59 The present invention preferably includes all functions in a miniature semiconductor chip, which as an integrated circuit, incorporates sensor, processing and transmitting units and control circuits.</p>

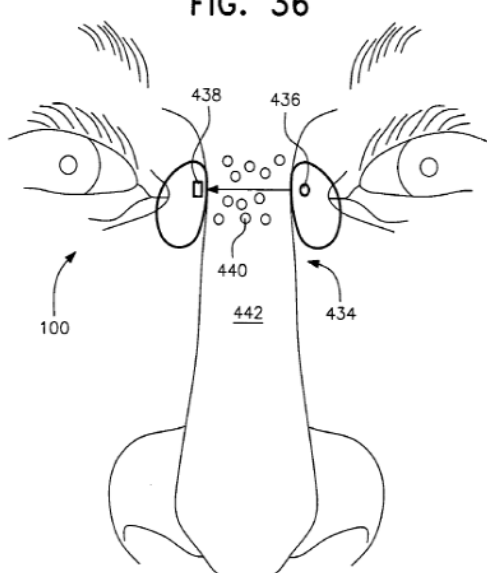
Claims	Prior Art
	<p data-bbox="583 180 890 207">Abreu '960 at 44:61–67</p> <p data-bbox="632 217 1974 358">The present invention includes means for directing preferably near-infrared energy into the surface of the skin at the end of the BTE tunnel, means for analyzing and converting the reflectance or back scattered spectrun into the concentration of the substance measured and support means for positioning the light source and detector means adjacent to the surface of the skin at the BTE tunnel site.</p> <p data-bbox="583 415 856 443">Irvin '332 at 2:20–34</p> <p data-bbox="632 453 1986 776">The wireless headset serves as a communication accessory to a base unit, such as a cellular telephone. When not in use, the wireless headset reverts to an inactive state and draws essentially no power. A power control circuit imparts automatic enabling capability to the wireless headset. Preferably, the power control circuit employs a sensor responsive to a natural stimulus associated with a person preparing to use the wireless headset. The wireless headset may sense motion, physical orientation, or user proximity. Sensing techniques include but are not limited to micro-electro-mechanical-systems (MEMS), mercury switches, attitude switches, and IR or acoustic proximity sensors. In some embodiments, the wireless headset transitions from the inactive state to an active state in response to a user picking it up or mounting the wireless headset upon their head.</p> <p data-bbox="583 833 840 860">Irvin '332 at 6:3–33</p> <p data-bbox="632 870 1986 1230">FIG. 4C illustrates the power control circuit 210 coupled with a contact switch functioning as the sensor 220. The contact switch includes a first and second contact connected to the power control circuit 210, and connecting member 222 actuated when the user places the wireless headset 100 on their body. When actuated, the connecting member 222 places the first and second contacts of sensor 220 in electrical connection. Ideally, the sensor 220, when configured as a contact sensor, is placed in the earpiece 130, in the retaining member 140, or in the retaining pad 150 shown in FIG. 1B. Preferably, when sensor 220 is configured as a contact-based proximity sensor, the contacting element of sensor 220 is sealed in a manner that minimizes contaminant ingress. In some embodiments, the contact sensor may be actuated by the user resiliently expanding the retaining member 140 or other element of the wireless headset 100 for mounting upon their body.</p> <p data-bbox="632 1240 1986 1451">FIG. 4D illustrates the power control circuit 210 coupled with a proximity sensor 220 not based on mechanical actuation. The proximity sensor 220 may be an active or passive IR sensor, may be a resistive or other type of thermal sensor, or may be an acoustic sensor. As a proximity sensor, sensor 220 provides indication that a user has placed the wireless headset 100 upon their body. As such, any of the foregoing or other presence sensing techniques may be used to implement sensor 220. Preferably, when configured as a proximity sensor, sensor 220 requires minimal power and provides a sensor output signal based on</p>

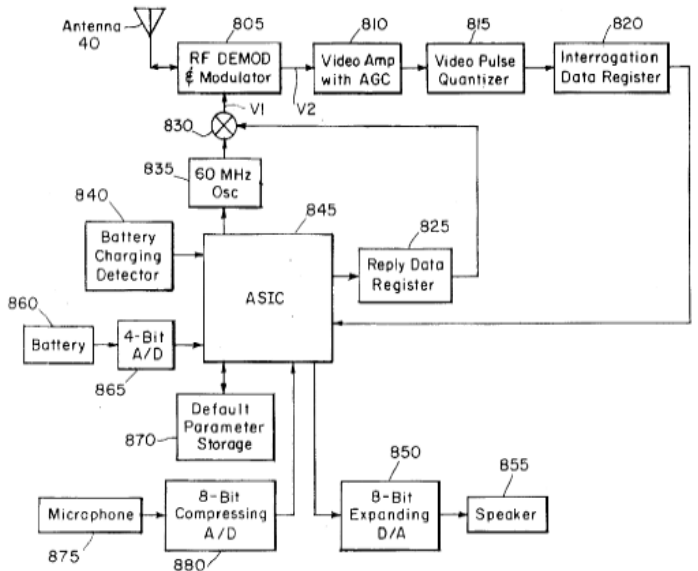
Claims	Prior Art
	<p>proximity to the user's body. Depending upon the sensing technology chosen, the control circuit 210 may be required to periodically stimulate sensor 220 to generate a detection response.</p> <p>Irvin '332 at 8:11–22</p> <p>FIGS. 6 and 7 illustrate operation of the wireless headset 100. FIG. 6 is a state diagram for one exemplary embodiment of the present invention. When the wireless headset 100 is not in use, it automatically enters an inactive state (610). As explained, this inactive state may constitute disconnecting the communications circuit 230 from the power source 240, or may simply represent placing the communications circuit 230 in a non-functional, low-power condition. Absent user stimulus, such as motion of the wireless headset 100 or placement of the wireless headset 100 upon the user's body, the wireless headset 100 remains in the inactive state (610).</p> <p>Luo '446 at [0034]</p> <p>In general, the location sensor module 18 and the hearing aid module 24 communicate via a bi-directional information signal 42. The hearing aid module 24 polls the location sensor module 18 via the bi-directional signal 42, on preferably a periodic basis, to determine whether the hearing aid 10 is in the ear of the hearing aid user. In this case, the bi-directional signal 42 acts as a polling signal. In response to polling done via the bi-directional information signal 42, the location sensor module 18 probes the outer environment of the hearing aid 10 and returns location information via the bi-directional location signal 42. In this case, the bi-directional signal 42 acts as a location information signal. If the location information indicates that the hearing aid 10 is in the ear of the hearing aid user (i.e. the in-the-ear case), the hearing aid 10 operates in full-function mode. If the location information indicates that the hearing aid 10 is not in the ear of the hearing aid user (i.e. the out-of-the-ear case), the hearing aid 10 operates in sleep mode. In sleep mode, the hearing aid 10 can essentially be considered to be off.</p> <p>Luo '446 at [0035]</p> <p>There are several general scenarios for the location of the hearing aid. In the first scenario, the hearing aid 10 is not in the ear of the user and the hearing aid 10 is in sleep mode. In this case, the location sensor module 18 is polled and the location information indicates that the hearing aid 10 is not in the ear of the hearing aid user and the hearing aid 10 continues to operate in sleep mode. In the second scenario, the hearing aid 10 has just been inserted into the ear of the hearing aid user and the hearing aid 10 was previously in sleep mode. In this case, the location sensor module 18 is polled and the location information indicates that the hearing aid 10 is in the ear of the hearing aid user. The hearing aid 10 then moves into full-function mode. In the third scenario, the hearing aid 10 is in full-function mode and is being taken out of the ear of the hearing aid user. The location sensor module 18 is polled and the location</p>

Claims	Prior Art
	<p>information indicates that the hearing aid 10 is no longer in the ear of the hearing aid user. The hearing aid 10 then moves into sleep mode.</p> <p>Luo '446 at [0036] The invention generally relies on shining light having a particular wavelength on human skin as well as the reflectance properties of human skin. The surface reflection of the human skin occurs at the surface of the epidermis and is approximately independent of the lighting wavelength and independent of human race. In the infrared (IR) range the wavelengths that can be used are greater than 800 nano-meters and the surface reflectance of human skin is more or less constant and close to 50%. The same hold true for orange and red light in the visible light range, namely wavelengths between approximately 600 and 800 nano-meters can be used, although the surface reflectance for darker skin is reduced for shorter wavelengths. Surface reflectance for dark skin is approximately 25% at 700 nano-meters. In the embodiments shown herein, the energy of choice in the sensor unit 18 is preferably infrared (IR) energy. However, long wavelength visible light energy can also be used as discussed below.</p> <p>Luo '446 at [0038] The optical window 56 is placed at a certain location on the shell 58 of the hearing aid 10. The location of the optical window 56 depends on whether the hearing aid is a BTE, ITE, ITC or CIC hearing aid. In addition, the location is chosen to minimize the distance between the optical window 56 and the skin 64. For optimal reflection of IR signals back through the optical window, the hearing aid user's skin (i.e. the reflecting surface) is required to be immediately over the optical window 56; otherwise the reflected IR energy will not be reflected back towards the reception unit 62. The optical window 56 is typically a small window having a diameter of approximately 1 mm for example. The optical window 56 can be made from IR grade glass or other suitable material that allows for the passage of infrared energy (a different material would be used if visible light is used rather than infrared light). The optical window 56 should be kept clean at all times in order to prevent the emitted infrared energy from being reflected back due to dirt and the like that may accumulate on the optical window 56 over time. The location of the optical window 56 will be discussed in further detail below.</p> <p>A POSITA would have been motivated to combine the teachings of the DHA system with Anderson and Abreu or Luo or Irvin to meet the claim limitations with a reasonable expectation of success for the reasons described above with respect to claim [1a] and [23a].</p> <p>To the extent that the DHA system, alone or in view of Anderson and Abreu or Luo or Irvin, does not disclose any elements of this claim limitation, it would have been obvious to a person having ordinary skill in the art</p>

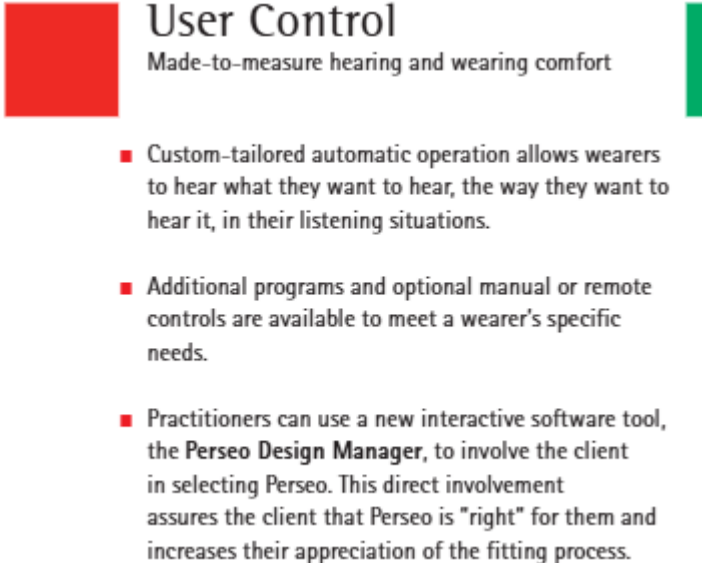
Claims	Prior Art
	<p>to combine, with the teachings of DHA systems, those elements, as known to them or as disclosed by other prior art identified in the cover pleading.</p>
<p>[56a] A personal audio system as recited in claim 42, wherein the first head-worn structure comprises: a radiation detector configured to capture radiation data; and</p>	<p>DHA systems, alone or in view of Anderson and Abreu or Luo or Irvin, expressly or inherently discloses “[a] personal audio system as recited in claim 42, wherein the first head-worn structure comprises: a radiation detector configured to capture radiation data.”</p> <p><i>See above at claim [23a].</i></p> <p>A POSITA would have been motivated to combine the teachings of the DHA system with Anderson and Abreu or Luo or Irvin to meet the claim limitations with a reasonable expectation of success for the reasons described above with respect to claim [1a] and [23a].</p> <p>To the extent that the DHA system, alone or in view of Anderson and Abreu or Luo or Irvin, does not disclose any elements of this claim limitation, it would have been obvious to a person having ordinary skill in the art to combine, with the teachings of DHA systems, those elements, as known to them or as disclosed by other prior art identified in the cover pleading.</p>
<p>[56b] first electronic circuitry provided within the first head-worn structure and electrically connected to the radiation detector, the first electronic circuitry configured to process at least the radiation data captured by the radiation detector.</p>	<p>DHA systems, alone or in view of Anderson and Abreu or Luo or Irvin, expressly or inherently discloses “first electronic circuitry provided within the first head-worn structure and electrically connected to the radiation detector, the first electronic circuitry configured to process at least the radiation data captured by the radiation detector.”</p> <p><i>See above at claim [23b].</i></p> <p>A POSITA would have been motivated to combine the teachings of the DHA system with Anderson and Abreu or Luo or Irvin to meet the claim limitations with a reasonable expectation of success for the reasons described above with respect to claim [1a] and [23a].</p> <p>To the extent that the DHA system, alone or in view of Anderson and Abreu or Luo or Irvin, does not disclose any elements of this claim limitation, it would have been obvious to a person having ordinary skill in the art to combine, with the teachings of DHA systems, those elements, as known to them or as disclosed by other prior art identified in the cover pleading.</p>
<p>57. A personal audio system as recited in claim 56, wherein the first electronic circuitry is configured to determine whether the first head-worn structure is</p>	<p>DHA systems, alone or in view of Anderson and Abreu or Luo or Irvin, expressly or inherently discloses “[a] personal audio system as recited in claim 56, wherein the first electronic circuitry is configured to determine whether the first head-worn structure is being worn by a user based on at least the processed radiation data.”</p> <p><i>See above at claim [23c].</i></p>

Claims	Prior Art
<p>being worn by a user based on at least the processed radiation data.</p>	<p>A POSITA would have been motivated to combine the teachings of the DHA system with Anderson and Abreu or Luo or Irvin to meet the claim limitations with a reasonable expectation of success for the reasons described above with respect to claim [1a] and [23a].</p> <p>To the extent that the DHA system, alone or in view of Anderson and Abreu, does not disclose any elements of this claim limitation, it would have been obvious to a person having ordinary skill in the art to combine, with the teachings of DHA systems, those elements, as known to them or as disclosed by other prior art identified in the cover pleading.</p>
<p>58. A personal audio system as recited in claim 57, wherein the radiation detector comprises a photodetector.</p>	<p>DHA systems, alone or in view of Abreu, expressly or inherently discloses “[a] personal audio system as recited in claim 57, wherein the radiation detector comprises a photodetector.”</p> <p>Abreu '960 at FIG. 35</p> <p>FIG. 35</p> 

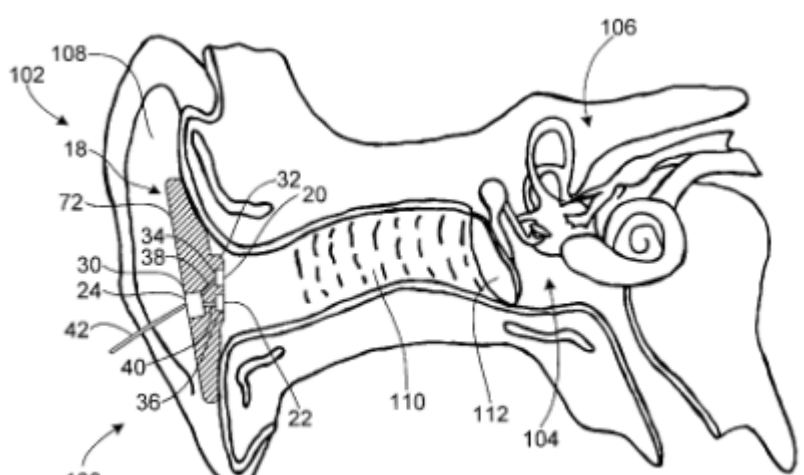
Claims	Prior Art
	<p>Abreu '960 at FIG. 36</p> <p style="text-align: center;">FIG. 36</p>  <p>Abreu '960 at 45:60–46:13</p> <p>In accordance with FIG. 35 a schematic view of one preferred reflectance measuring apparatus of the present invention is shown. FIG. 35 shows a light source 420 such as an infrared LED and a photodetector 422 located side-by-side and disposed within support structure 426 such as a medial canthal pad or modified nose pads of eyeglasses directing radiation 424 at the BTE tunnel 430 with said light source 420 laying in apposition to the skin 428 at the BTE tunnel 430. The light source 420 delivers the radiation 424 to the skin 428 at the BTE tunnel which is partially absorbed according to the interaction with the substance 432 being measured resulting in attenuated radiation 425. Part of the radiation 424 is then absorbed by the substance 432 and the resulting radiation 425 emitted from BTE tunnel 430 is collected by the photodetector 422 and converted by a processor into the blood concentration of the substance 432. Thin skin 428 is the only tissue interposed between radiation 424, 425 and the substance 432 being measured. The concentration of the substance 432 is accomplished by detecting the magnitude of light attenuation collected which is caused by the absorption signature of the substance being measured.</p>

Claims	Prior Art
	<p>A POSITA would have been motivated to combine the teachings of the DHA system with Abreu to meet the claim limitations with a reasonable expectation of success for the reasons described above with respect to claim [1a] and [23a].</p> <p>To the extent that the DHA system, alone or in view of Abreu, does not disclose any elements of this claim limitation, it would have been obvious to a person having ordinary skill in the art to combine, with the teachings of DHA systems, those elements, as known to them or as disclosed by other prior art identified in the cover pleading.</p>
<p>59. A personal audio system as recited in claim 42, wherein the first head-worn structure comprises first electronic circuitry configured to provide audio alerts to the user via the first speaker.</p>	<p>DHA systems, alone or in view of Anderson, expressly or inherently discloses “[a] personal audio system as recited in claim 42, wherein the first head-worn structure comprises first electronic circuitry configured to provide audio alerts to the user via the first speaker.”</p> <p>Phonak – Fitting Guide for MLxi and MLxi Baha at Page 5</p> <p>Beep type MLxi/MLxi Baha plays a confirmation beep at start up or when the system synchronizes with a transmitter. You can change the type of beep tone using the drop-down menu.</p> <p>Anderson '783 at FIG. 8</p>  <p>The diagram (FIG. 8) illustrates a personal audio system. It features an Antenna (40) connected to an RF DEMOD & Modulator (805). The RF DEMOD & Modulator (805) is connected to a Video Amp with AGC (810), which in turn connects to a Video Pulse Quantizer (815). The Video Pulse Quantizer (815) outputs to an Interrogation Data Register (820). A 60 MHz Osc (835) is connected to a multiplier (830) and the ASIC (845). The multiplier (830) also receives input from the Interrogation Data Register (820) via a signal labeled V2. The ASIC (845) is a central component connected to several other blocks: a Battery Charging Detector (840), a 4-Bit A/D (860), a Default Parameter Storage (870), an 8-Bit Compressing A/D (880), an 8-Bit Expanding D/A (850), and a Reply Data Register (825). The 8-Bit Expanding D/A (850) is connected to a Speaker (855). A Battery (865) is connected to the 4-Bit A/D (860).</p>

Claims	Prior Art
	<p>Anderson '783 at 3:40–41 FIG. 8 is a block diagram showing details of an earpiece for a full-featured preferred embodiment of the invention.</p> <p>Anderson '783 at 24:15 Unsolicited information:</p> <p>Anderson '783 at 24:15–26 The invention provides useful unsolicited information (i.e., information automatically provided to the user by the system, not as a result of a user request for information) to the user of the system in a manner that can be made imperceptible to a casual observer if desired. Information about the status of the hearing aid system itself, such as a verbal warning that an earpiece or RPU battery voltage is low, is provided to the user as detailed earlier. Automatic measurements of heart pulse rate, body temperature and walking distance traveled can be taken by sensors 950 in the RPU 16 or earpiece 10 and stored in the RPU DSP 948.</p> <p>A POSITA would have been motivated to combine the teachings of the Phonak system with Anderson '783 to meet the claim limitations with a reasonable expectation of success for the reasons described above with respect to claim 6.</p> <p>To the extent that the Phonak system and Anderson '783, a does not disclose any elements of this claim limitation, it would have been obvious to a person having ordinary skill in the art to combine, with the teachings of the Phonak system and Anderson '783 or Jannard, those elements, as known to them or as disclosed by other prior art identified in the cover pleading.</p>
<p>60. A personal audio system as recited in claim 59, wherein the audio alerts are dependent on one or more user preference settings.</p>	<p>DHA systems, alone or in view of Anderson, expressly or inherently discloses “[a] personal audio system as recited in claim 59, wherein the audio alerts are dependent on one or more user preference settings.”</p>

Claims	Prior Art
	<p>Phonak – Perseo Brochure at Page 7</p>  <p>User Control Made-to-measure hearing and wearing comfort</p> <ul style="list-style-type: none"> ■ Custom-tailored automatic operation allows wearers to hear what they want to hear, the way they want to hear it, in their listening situations. ■ Additional programs and optional manual or remote controls are available to meet a wearer's specific needs. ■ Practitioners can use a new interactive software tool, the Perseo Design Manager, to involve the client in selecting Perseo. This direct involvement assures the client that Perseo is "right" for them and increases their appreciation of the fitting process. <p>Anderson '783 at 24:15 Unsolicited information:</p> <p>Anderson '783 at 24:15–26 The invention provides useful unsolicited information (i.e., information automatically provided to the user by the system, not as a result of a user request for information) to the user of the system in a manner that can be made imperceptible to a casual observer if desired. Information about the status of the hearing aid system itself, such as a verbal warning that an earpiece or RPU battery voltage is low, is provided to the user as detailed earlier. Automatic measurements of heart pulse rate, body temperature and walking distance traveled can be taken by sensors 950 in the RPU 16 or earpiece 10 and stored in the RPU DSP 948.</p> <p>A POSITA would have been motivated to combine the teachings of the Phonak system with Anderson '783 to meet the claim limitations with a reasonable expectation of success for the reasons described above with respect to claim 6.</p>

Claims	Prior Art
	<p>To the extent that the Phonak system and Anderson '783 or Jannard, a does not disclose any elements of this claim limitation, it would have been obvious to a person having ordinary skill in the art to combine, with the teachings of the Phonak system and Anderson '783 or Jannard, those elements, as known to them or as disclosed by other prior art identified in the cover pleading.</p>
<p>[61a] A personal audio system as recited in claim 42, wherein the first head-worn structure of the head-worn wearable electronic device comprises: a first microphone configured to pick up at least audio input from at least the user, the first microphone being at least partially in the first head-worn structure; and</p>	<p>DHA systems, alone or in view of Anderson and Nemirovski, expressly or inherently discloses “[a] personal audio system as recited in claim 42, wherein the first head-worn structure of the head-worn wearable electronic device comprises: a first microphone configured to pick up at least audio input from at least the user, the first microphone being at least partially in the first head-worn structure.”</p> <p><i>See above at claim 6 and [50a].</i></p> <p>A POSITA would have been motivated to combine the teachings of the DHA system with Anderson and Nemirovski to meet the claim limitations with a reasonable expectation of success for the reasons described above with respect to claims [1a] and 11.</p> <p>To the extent that the DHA system, alone or in view of Anderson and Nemirovski, does not disclose any elements of this claim limitation, it would have been obvious to a person having ordinary skill in the art to combine, with the teachings of DHA systems, those elements, as known to them or as disclosed by other prior art identified in the cover pleading.</p>
<p>[61b] a second microphone configured to pick up at least ambient sound input associated with an environment in the vicinity of the first head-worn structure,</p>	<p>DHA systems, alone or in view of Anderson and Nemirovski, expressly or inherently discloses “a second microphone configured to pick up at least ambient sound input associated with an environment in the vicinity of the first head-worn structure.”</p> <p><i>See above at claim 6 and [50b].</i></p> <p>A POSITA would have been motivated to combine the teachings of the DHA system with Anderson and Nemirovski to meet the claim limitations with a reasonable expectation of success for the reasons described above with respect to claims [1a] and 11.</p> <p>To the extent that the DHA system, alone or in view of Anderson and Nemirovski, does not disclose any elements of this claim limitation, it would have been obvious to a person having ordinary skill in the art to combine, with the teachings of DHA systems, those elements, as known to them or as disclosed by other prior art identified in the cover pleading.</p>

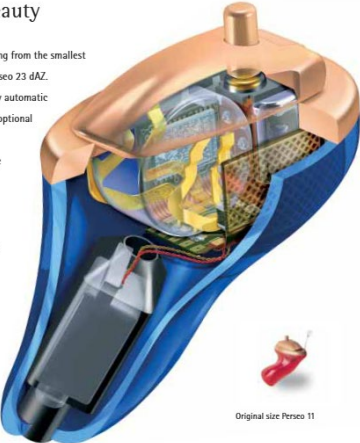
Claims	Prior Art
<p>[61c] wherein the second microphone is positioned to be further away from mouth of the user than the first microphone when at least the first head-worn structure of the head-worn wearable electronic device is being worn.</p>	<p>DHA systems, alone or in view of Nemirovski, expressly or inherently discloses “wherein the second microphone is positioned to be further away from mouth of the user than the first microphone when at least the first head-worn structure of the head-worn wearable electronic device is being worn.”</p> <p>Nemirovski '379 at FIG. 3</p>  <p style="text-align: center;">FIG. 3</p> <p>Nemirovski '379 at 2:54–56</p> <p>FIG. 3 is a cross-sectional view of an ear microphone assembly, or earset, disposed with respect to the ear of a person, and the ear also shown in cross-section;</p> <p>Nemirovski '379 at 6:36–41</p> <p>The microphone 20 is disposed on an interior portion 32 of the housing 30. Accordingly, the interior portion 32 is sized and shaped to position the microphone 20 with respect to the ear 100 and/or ear canal 110 as desired. In one embodiment, the interior portion 32 positions the microphone 20 and/or the speaker 22 in the concha 122.</p> <p>Nemirovski '379 at 9:19–25</p> <p>With additional reference to FIG. 5, in another embodiment of the present invention, the earset 18 can include an external microphone 70 in addition to the microphone 20 directed toward the ear canal of the user (also referred to herein as an internal microphone 20). The external microphone 70 detects air</p>

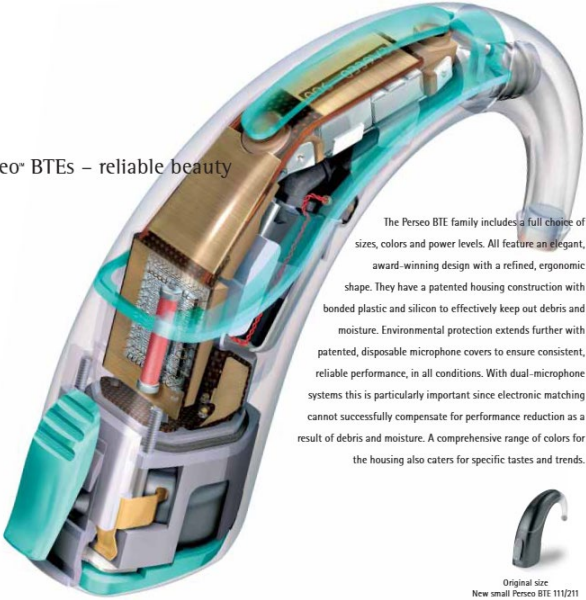
Claims	Prior Art
	<p>pressure changes, or sounds, originating outside the user's ear.</p> <p>Nemirovski '379 at 9:33–39 The external microphone 70 is disposed on a exterior portion 72 of the housing 30 and is pointed away from the wearer of the earset 18 assembly. The external microphone 70 can be secured to the housing 30 with adhesive and can be disposed in a recess defined by the exterior portion 72. Similar to the internal microphone 20, the external microphone 70 is small (e.g., less than two mm in size).</p> <p>A POSITA would have been motivated to combine the teachings of the DHA system with Nemirovski to meet the claim limitations with a reasonable expectation of success for the reasons described above with respect to claims [1a] and 11.</p> <p>To the extent that the DHA system, alone or in view of Nemirovski, does not disclose any elements of this claim limitation, it would have been obvious to a person having ordinary skill in the art to combine, with the teachings of DHA systems, those elements, as known to them or as disclosed by other prior art identified in the cover pleading.</p>
<p>[62] A personal audio system as recited in claim 42, wherein the first head-worn structure comprises: a physical condition sensor configured to acquire data pertaining to a physical condition of the user, the physical condition sensor being provided at least partially internal to the first head-worn structure.</p>	<p>DHA systems, alone or in view of Anderson, expressly or inherently discloses “[a] personal audio system as recited in claim 42, wherein the first head-worn structure comprises: a physical condition sensor configured to acquire data pertaining to a physical condition of the user, the physical condition sensor being provided at least partially internal to the first head-worn structure.”</p> <p>Anderson '783 at 2:40–59 Another object of this invention is to provide supplemental audio information (e.g. a verbal warning from the RPU that an earpiece battery is low) and communication services (e.g., cellular telephone and paging services) to the user via the hearing aid system. Such services are accessed in a manner that can be made imperceptible to a casual observer if desired. The RPU (which contains a digital signal processor or other computer) acts as one source of information, e.g., by using a synthesized voice message to provide the time of day, and the RPU can also be used to access a secondary wireless link to the general subscriber telephone network or voice paging services. User control of hearing aid parameters and requests for information are accomplished using pushbuttons located on the RPU (including pushbuttons suitable for data entry in a covert manner) or voice recognition. Many sensors and peripheral devices can reside in or be attached to the RPU by wired or wireless means, and can provide a variety of information for different applications (e.g., heart pulse rate) as audio in the user's earpiece.</p>

Claims	Prior Art
	<p>Anderson '783 at 24:15–26</p> <p>The invention provides useful unsolicited information (i.e., information automatically provided to the user by the system, not as a result of a user request for information) to the user of the system in a manner that can be made imperceptible to a casual observer if desired. Information about the status of the hearing aid system itself, such as a verbal warning that an earpiece or RPU battery voltage is low, is provided to the user as detailed earlier. Automatic measurements of heart pulse rate, body temperature and walking distance traveled can be taken by sensors 950 in the RPU 16 or earpiece 10 and stored in the RPU DSP 948.</p> <p>Anderson '783 at 25:43–50</p> <p>In emergency situations, the system may automatically dial an appropriate emergency service (based on information previously stored by the user in the RPU DSP 948), send voice or digital information regarding location, current medical situation (e.g. heart pulse rate) and medical history data (stored in the RPU DSP 948) in addition to providing normal telephone voice communication via the secondary wireless link circuitry 944.</p> <p>To the extent that the DHA system, alone or in view of Anderson, does not disclose any elements of this claim limitation, it would have been obvious to a person having ordinary skill in the art to combine, with the teachings of DHA systems, those elements, as known to them or as disclosed by other prior art identified in the cover pleading.</p>
<p>63. A personal audio system as recited in claim 62, wherein the acquired data by the physical condition sensor includes data pertaining to or dependent on heart beats of the user.</p>	<p>DHA systems, alone or in view of Anderson, expressly or inherently discloses “[a] personal audio system as recited in claim 62, wherein the acquired data by the physical condition sensor includes data pertaining to or dependent on heart beats of the user.”</p> <p>Anderson '783 at 2:40–59</p> <p>Another object of this invention is to provide supplemental audio information (e.g. a verbal warning from the RPU that an earpiece battery is low) and communication services (e.g., cellular telephone and paging services) to the user via the hearing aid system. Such services are accessed in a manner that can be made imperceptible to a casual observer if desired. The RPU (which contains a digital signal processor or other computer) acts as one source of information, e.g., by using a synthesized voice message to provide the time of day, and the RPU can also be used to access a secondary wireless link to the general subscriber telephone network or voice paging services. User control of hearing aid parameters and requests for information are accomplished using pushbuttons located on the RPU (including pushbuttons suitable for data entry in a covert manner) or voice recognition. Many sensors and peripheral devices can reside in or</p>

Claims	Prior Art
	<p>be attached to the RPU by wired or wireless means, and can provide a variety of information for different applications (e.g., heart pulse rate) as audio in the user's earpiece.</p> <p>Anderson '783 at 24:15–26</p> <p>The invention provides useful unsolicited information (i.e., information automatically provided to the user by the system, not as a result of a user request for information) to the user of the system in a manner that can be made imperceptible to a casual observer if desired. Information about the status of the hearing aid system itself, such as a verbal warning that an earpiece or RPU battery voltage is low, is provided to the user as detailed earlier. Automatic measurements of heart pulse rate, body temperature and walking distance traveled can be taken by sensors 950 in the RPU 16 or earpiece 10 and stored in the RPU DSP 948.</p> <p>To the extent that the DHA system, alone or in view of Anderson, does not disclose any elements of this claim limitation, it would have been obvious to a person having ordinary skill in the art to combine, with the teachings of DHA systems, those elements, as known to them or as disclosed by other prior art identified in the cover pleading.</p>
<p>64[pre]. A personal audio system as recited in claim 24, wherein the head-worn wearable electronic device comprises:</p>	<p>DHA systems, alone or in view of Anderson, expressly or inherently discloses “[a] personal audio system as recited in claim 24, wherein the head-worn wearable electronic device comprises.”</p> <p><i>See</i> above at claim 1[pre]; 24.</p> <p>A POSITA would have been motivated to combine the teachings of the DHA system with Anderson and Nemirovski to meet the claim limitations with a reasonable expectation of success for the reasons described above with respect to claims [1a] and 11.</p> <p>To the extent that the DHA system, alone or in view of Anderson and Nemirovski, does not disclose any elements of this claim limitation, it would have been obvious to a person having ordinary skill in the art to combine, with the teachings of DHA systems, those elements, as known to them or as disclosed by other prior art identified in the cover pleading.</p>
<p>[64a] a head-worn structure including at least a first microphone, a second microphone, a speaker and a processing unit;</p>	<p>DHA systems, alone or in view of Anderson, expressly or inherently discloses “a head-worn structure including at least a first microphone, a second microphone, a speaker and a processing unit.”</p>

Claims	Prior Art
	<p data-bbox="583 175 1150 212">Phonak - SmartLink SX User Guide at p. 32</p>  <p data-bbox="583 672 1010 709">Phonak – Perseo Brochure at p. 8</p> 

Claims	Prior Art
	<p>Phonak – Perseo Brochure at Page 22</p> <h1 data-bbox="655 224 1255 373">The beauty of high technology</h1> <p data-bbox="703 378 957 402">Perseo[®] ITEs – hidden beauty</p> <div data-bbox="630 487 655 535">2</div> <p data-bbox="703 427 1003 589">The Perseo family includes a full line of ITEs ranging from the smallest Perseo 11 CIC to the elegant dual-microphone Perseo 23 dAZ. All CIC and ITE instruments, like the BTEs, are fully automatic and have PersonalLogic. A TacTronic switch is optional on all ITE models and easily added after fitting by exchanging battery covers. This is a very simple and practical means of providing extra control and additional programs.</p> <p data-bbox="703 617 915 695">In many countries, ITEs are available with Phonak eShell. These electronically created shells are hypoallergenic, exceptionally durable and more easily retained in the wearer's ear.</p>  <p data-bbox="1150 787 1234 799">Original size Perseo 11</p>

Claims	Prior Art
	<p data-bbox="583 175 1062 207">Phonak – Perseo Brochure at Page 23</p> <div data-bbox="627 224 1255 824">  <p data-bbox="632 381 905 406">Perseo® BTEs – reliable beauty</p> <p data-bbox="968 440 1255 691">The Perseo BTE family includes a full choice of sizes, colors and power levels. All feature an elegant, award-winning design with a refined, ergonomic shape. They have a patented housing construction with bonded plastic and silicon to effectively keep out debris and moisture. Environmental protection extends further with patented, disposable microphone covers to ensure consistent, reliable performance, in all conditions. With dual-microphone systems this is particularly important since electronic matching cannot successfully compensate for performance reduction as a result of debris and moisture. A comprehensive range of colors for the housing also caters for specific tastes and trends.</p> <p data-bbox="1136 800 1247 824">Original size New small Perseo BTE 111/211</p> </div>

Claims

Prior Art

Phonak – Perseo Brochure at Page 26

Perseo 11

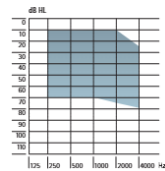


||| Digital microphone

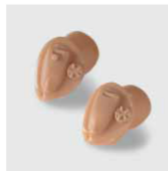
OPTIONS
||| Tactronic switch
||| Control beeps

Ear Simulator data
45 dB max. gain
119 dB SPL max. output

2 cc data
35 dB max. gain
108 dB SPL max. output



Perseo 12

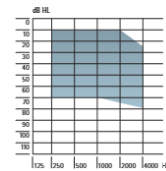


||| Digital microphone
||| Volume control

OPTIONS
||| Tactronic switch
||| Control beeps

Ear Simulator data
50 dB max. gain
119 dB SPL max. output

2 cc data
40 dB max. gain
108 dB SPL max. output



Perseo 22

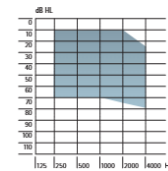


||| Telecoil

OPTIONS
||| Tactronic switch
||| Control beeps
||| Remote control WatchPlot2
and SoundPlot

Ear Simulator data
58 dB max. gain
118 dB SPL max. output

2 cc data
48 dB max. gain
108 dB SPL max. output



Claims

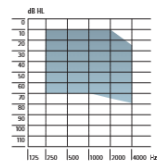
Prior Art

Phonak – Perseo Brochure at Page 27

Perseo 23 dAZ



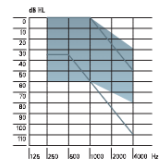
- Adaptive digital AudioZoom
- Telecoil
- OPTIONS
 - Tactronic switch
 - Control beeps
 - Remote control WatchPlot2 and SoundPilot
- Ear Simulator data
 - 58 dB max. gain
 - 118 dB SPL max. output
- 2 cc data
 - 45 dB max. gain
 - 108 dB SPL max. output



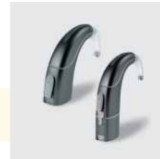
Perseo 111 dAZ
Perseo 111 dAZ FM



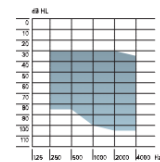
- FM solutions
- Adaptive digital AudioZoom
- Telecoil
- Tactronic switch
- OPTIONS
 - Control beeps
 - Remote control WatchPlot2 and SoundPilot
- Ear Simulator data
 - 54 dB max. gain
 - 127 dB SPL max. output
- 2 cc data
 - 47 dB max. gain
 - 121 dB SPL max. output



Perseo 211 dAZ
Perseo 211 dAZ FM



- FM solutions
- Adaptive digital AudioZoom
- Telecoil
- Tactronic switch
- OPTIONS
 - Control beeps
 - Remote control WatchPlot2 and SoundPilot
- Ear Simulator data
 - 62 dB max. gain
 - 130 dB SPL max. output
- 2 cc data
 - 55 dB max. gain
 - 125 dB SPL max. output



Technical Data

Claims

Prior Art

Phonak – Perseo Brochure at Page 28

Perseo 311 dAZ Forte
Perseo 311 dAZ Forte FM



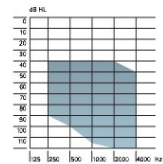
- FM solutions
- Adaptive digital AudioZoom
- Telcoil
- Tachonic Switch

OPTIONS

- Control beeps
- Remote control WatchPilot2 and SoundPilot

Ear Simulator data
72 dB max. gain
135 dB SPL max. output

2 cc data*
65 dB max. gain
130 dB SPL max. output



Remote controls



WatchPilot™2



SoundPilot™



FM Transmitters



HandyMic™



TelCom™



Campus S™

Valente - Effect of type of noise and loudspeaker array on the performance at 398

Differences in performance between omnidirectional and directional microphones were evaluated between two loudspeaker conditions (single loudspeaker at 180°; diffuse using eight loudspeakers set 45° apart) and two types of noise (steady-state HINT noise; R-Space™ restaurant noise). Twenty-five participants were fit bilaterally with Phonak Perseo hearing aids using the manufacturer's recommended procedure. After wearing the hearing aids for one week, the

Valente - Effect of type of noise and loudspeaker array on the performance at 399

Recently, Phonak introduced its digital hearing aid in which one feature is an adaptive directional microphone. This adaptive microphone reportedly ensures identical sensitivity and phase characteristics of the dual microphones and a mechanism to *automatically* change the polar pattern from omnidirectional to any directional pattern while compensating for the changes in low-frequency response associated with each pattern. This hearing aid can be programmed

Claims	Prior Art
	<p data-bbox="583 175 1514 212">Phonak – The Hearing Review Article re Perseo (Sep. 5, 2003) at Page 2</p> <p data-bbox="636 228 1226 521">Perseo is a 20-channel digital hearing instrument. The default base-program is optimized for speech understanding in quiet environments. A second base-program with different frequency response settings and optional adaptive directional microphones provides additional support in noisy situations. The acoustic environment is continuously monitored and the appropriate base-program is automatically selected. This dual-base automatic function allows for independent fine tuning in one listening situation without compromising performance in the other. The automatic performance is designed to be easily personalized for the individual needs and preferences of the user.</p> <p data-bbox="583 545 1176 583">Phonak – Annual Report (2003–2004) at p. 12</p> <p data-bbox="636 599 1087 989">Within the FM business segment, the Group continued to leverage existing key technology into new applications. With SmartLink, Phonak has achieved another industry first by allowing hearing instrument users to access other communications systems. SmartLink combines a remote control, a transmitter, a directional microphone system and a Bluetooth link to mobile phone all in one single device. Phonak also presented EduLink, an innovative FM device to assist the learning development of children with Auditory Processing Disorder (APD). EduLink was perceived very positively within the medical community and recent classroom tests show favourable results.</p>

Claims

Prior Art

Phonak – Collections 2002 Product List at Page 6

Elite



Digital Perception Processing (DPP)



AutoSelect (AS)

New Phonak sound:
unmatched clarity
Digital Perception Processing (DPP) is based on a cochlear sound analysis across 20 critical bands – exactly like the analysis performed by the healthy human ear. But Clari gets a step further than previous overlapping multi-band filter systems. DPP models the active, load-ress-dependent cochlear signal processing in all 20 critical bands.

20 critical bands:
overlapping and dynamic
The cochlear filters are overlapping, but the degree of overlap is dependent on the frequency, level, and bandwidth of the sound signal. This is precisely how Clari's filters work. Digital Perception Processing computes the cochlear excitation patterns and thus determines sound perception. DPP is the first to mimic the healthy human ear in this way.

Automatic, intelligent program selection
Clari listens and automatically activates the optimal DPP algorithm for the acoustic environment: Program 1 with QuietAdapt for quiet situations or Program 2 with NoiseAdapt and noise control for noisy environments.

In quiet situations
QuietAdapt DPP is designed for clear hearing in quiet surroundings. QuietAdapt controls gain in 20 critical bands to correct for loudness and provide a natural sound quality. By preserving the perception patterns of complex sounds, QuietAdapt assures that low-level signal components are not over-

amplified to the point where they would adversely affect speech intelligibility. This allows Clari to deliver only those sound components that would be naturally audible.

In noise
NoiseAdapt DPP is the ideal program for maximizing audibility in noisy environments. The special NoiseAdapt algorithm works with the Fine-scale Noise Canceller – and if available with the Adaptive digital AudioZoom feature – to optimize speech understanding in noise.



Claims

Prior Art

Phonak – Collections 2002 Product List at Page 7



Adaptive digital AudioZoom (dAZ)



Fine-scale Noise Canceler (FNC)

Active and adaptive separation of speech from noise
The Speech-to-Noise Ratio (SNR) can only be optimised if the noise is spatially separated from the speech signal with directional microphone technology. As a pioneer in multi-microphone technologies, Phonak introduced AudioZoom, a feature which significantly improves listening and understanding in noisy situations, enhancing user satisfaction. Now, Phonak goes a step further in the refinement of multi-microphone technology with Adaptive digital AudioZoom: it seeks, finds, and suppresses noise automatically.

Fully automatic, fully adaptive
With its unique spatial processing system, Phonak's Adaptive digital AudioZoom can detect noise not only from directly behind, but anywhere in the near plane. Whether the hearing-impaired person is moving or the location of the noise source changes, the directional characteristic is constantly adapted to the listening situation. In this manner, the Adaptive digital AudioZoom automatically establishes the best possible SNR.

Twenty-band noise control
Claro constantly monitors the Speech-to-Noise Ratio (SNR) in all 20 bands. If a poor SNR is measured in certain bands, the Fine-scale Noise Canceler will reduce gain only in those bands. The result of this precision noise control is a selective and intelligent improvement of the SNR with only a marginal influence on the speech signal.

Intelligent
The magnitude of gain reduction is not constant. In frequency bands with higher SNRs and in bands which are important for speech intelligibility, gain is reduced only slightly. This intelligent control algorithm reduces the masking and distractive effects of background noise. Important speech information is preserved.

Precise and selective
All frequency bands are separately analysed and controlled. Gain is reduced selectively only in those bands where noise is present. This high-resolution analysis and gain control makes it possible to suppress noise with extremely high accuracy. Neighboring frequency bands without noise remain unaffected.

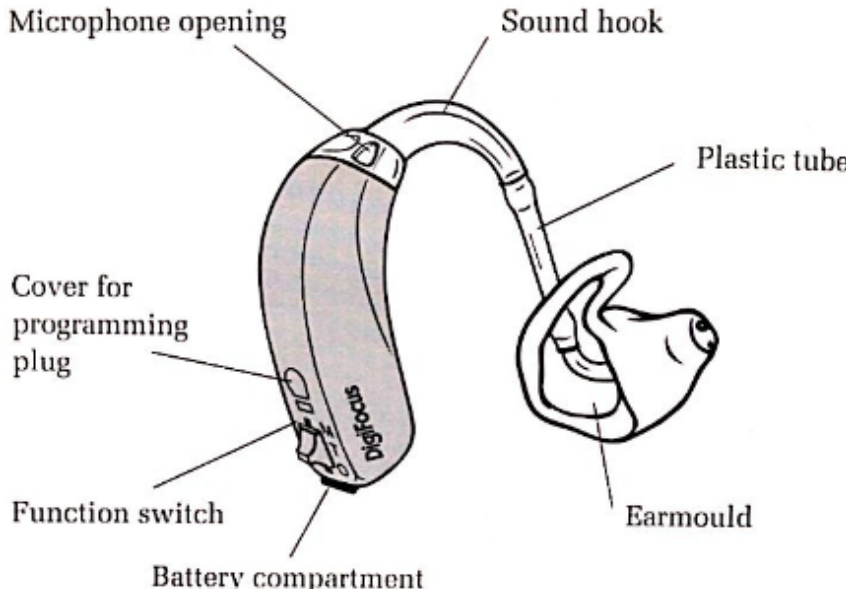
Claims

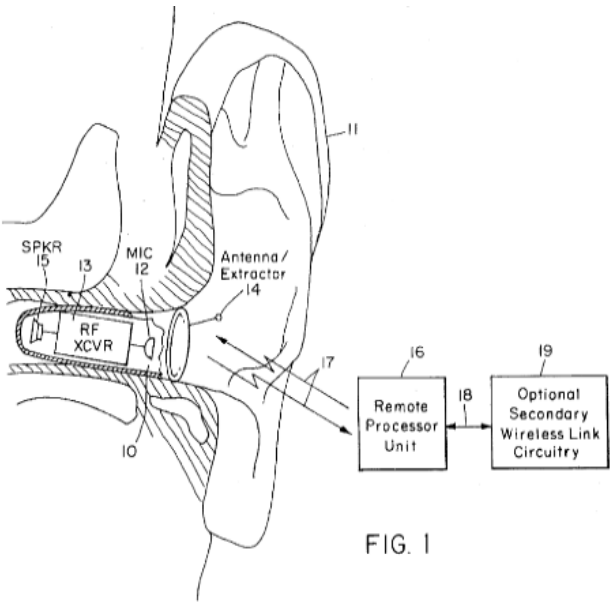
Prior Art

Phonak – Collections 2002 Product List at Page 8

	Claro 11 (CIC)	Claro 21 dAZ	Claro 22
ITE instruments			
			
Peak values G max./OSPL 90	45 dB/117 dB SPL	58 dB/117 dB SPL	58 dB/117 dB SPL
Features	World's first digital microphone TacTronic switch standard	In-the-Ear hearing computer with Adaptive digital Audio/oom (dAZ)	Semi-modular Claro In-the-Canal instrument
Description	Digital Perception Processing (DPP) AutoSelect (AS) Fine-scale Noise Canceler (FNC) 3 hearing programs + "mute" Battery size 10 TacTronic program selection switch	Digital Perception Processing (DPP) AutoSelect (AS) Fine-scale Noise Canceler (FNC) Adaptive digital Audio/oom (dAZ) 3 hearing programs + "mute" Battery size 312 Telecoil	Digital Perception Processing (DPP) AutoSelect (AS) Fine-scale Noise Canceler (FNC) 3 hearing programs + "mute" Battery size 312 Telecoil
Options	Rotary volume control potentiometer Colors: Pink, tan, cocoa, medium brown	Remote controls Watch/Plot and Sound/Plot FM TacTronic program selection switch On/Off switch (battery cover) Colors: Pink, tan, cocoa, medium brown	Remote controls Watch/Plot and Sound/Plot FM TacTronic program selection switch On/Off switch (battery cover) Colors: Pink, tan, cocoa, medium brown
Programming	Loudness Perception Profile (LPP) with Active Response Unit (ARU) Fitting formulas: Phonak Digital, Phonak Digital Ski Slope and DSL (i/o)	Loudness Perception Profile (LPP) with Active Response Unit (ARU) Fitting formulas: Phonak Digital, Phonak Digital Ski Slope and DSL (i/o)	Loudness Perception Profile (LPP) with Active Response Unit (ARU) Fitting formulas: Phonak Digital, Phonak Digital Ski Slope and DSL (i/o)

Claims	Prior Art
	<p>Lewis – Speech perception in noise: directional microphones (2004) at Page 431</p> <p>All subjects were fit with digital Phonak Claro 311 dAZ BTE hearing aids bilaterally. All earmolds had select-a-vent (SAV) venting and #13 or 3 mm horn tubing. In addition to the hearing aids, subjects were fit with Phonak Microlink ML8 FM receivers bilaterally. These FM receivers attach to the bottom of a BTE hearing aid and may be utilized in either the “FM only” mode, which attenuates the hearing aid microphone by 20 dB, or in the FM plus hearing aid mode, which allows for FM input and input of environmental sounds via the hearing aid microphones simultaneously at the same output level. The Phonak TX3 HandyMic FM transmitter served as the FM transmitter. This transmitter has three microphone options: (1) Wide Angle, which picks up sounds arriving from all directions around the transmitter microphone equally; (2) Zoom, which provides reduced amplification to signals arriving from the rear (cardioid); and (3) SuperZoom, which provides reduced amplification to signals originating from the rear and the sides (hypercardioid). The</p>

Claims	Prior Art
	<p data-bbox="583 175 1449 207">Oticon – DigiFocus Compact User Instructional Brochure at Page 4</p>  <p data-bbox="583 844 1075 876">Oticon – DigiFocus Description at p. 1</p> <p data-bbox="634 896 1312 966">DigiFocus is available both as behind-the-ear models (BTE) in colours which match the hair tone, and as in-the-ear (ITE) models which fit discreetly into the ear.</p>  <p data-bbox="714 1291 1239 1315"><i>DigiFocus Compact, Canal, Half Shell and Low Profile.</i></p>

Claims	Prior Art
	<p data-bbox="583 175 1220 207">Oticon – DigiFocus White Paper (1996) at Page 2</p> <p data-bbox="636 224 1310 358">Manufacturers of hearing instruments have worked hard to overcome the troubles of size and battery consumption. With the thickness of digital microchips at only 1/100 of a human hair, this has now been accomplished in the world's first fully digital, at the ear hearing instrument - DigiFocus. Inside DigiFocus there is computing power available which equals that of a PC.</p> <p data-bbox="583 386 1234 418">Oticon – DigiFocus White Paper (1996) at Page 11</p> <p data-bbox="636 435 1310 570">Many people find it uncomfortable to be wired to a computer. This "constricted" feeling can be avoided by using the OtiSet programming belt, which is worn on the shoulder. The programming belt provides cordless infrared communication between the PC and the hearing instrument while in the ear. Alternatively, you may use the HI-PRO interface for a wired solution.</p> <p data-bbox="636 591 1310 748">OtiSet offers a choice of three different views of the fitting. The Audiological View provides a visual guide to instrument performance in four typical sound environments. The Technical View illustrates traditional technical information such as coupler data and insertion gain at different input levels. With the Binaural View you are able to examine settings for both right and left instruments at the same time.</p> <p data-bbox="583 776 894 808">Anderson '783 at FIG. 1</p>  <p data-bbox="995 1349 1062 1382">FIG. 1</p>

Claims	Prior Art
	<p data-bbox="583 180 919 207">Anderson '783 at 3:13–18</p> <p data-bbox="632 217 1969 358">FIG. 1 is a block diagram system overview illustrating how an earpiece (comprising a microphone, RF transceiver, dual-use antenna/extractor and speaker) is worn in the ear, communicates via a primary two-way RF link with an RPU, and also communicates via an optional secondary wireless link to a telephone system.</p> <p data-bbox="583 415 945 443">Anderson '783 at 3:52–4:14</p> <p data-bbox="632 453 1986 1068">One preferred embodiment of the invention herein described is shown in block diagram form in FIG. 1. In this figure, the earpiece 10 shown worn in the ear 11 uses a standard completely-in-the-canal (CIC) housing, well known in the art, although many other earpiece housing types known in the art (e.g., behind-the-ear, or BTE) can be used. The earpiece shown in FIG. 1 comprises a microphone 12, RF transceiver 13 with an antenna that doubles as an earpiece extractor 14, and speaker 15. Although the dual-use antenna/extractor arrangement shown in FIG. 1 is often desirable, a separate antenna or antennae can be used in conjunction with many different types of earpiece extractors as appropriate for a particular application. Antennae and extractors can be disguised as ear hair or jewelry (e.g., earrings), or conductive filaments can be permanently implanted in the ear cartilage to act as antennae using a process similar to ear piercing. The earpiece 10 communicates with a remote processor unit (RPU) 16 via a primary two-way RF link 17, although many other wireless link media (e.g., ultrasonic or infrared) may be used instead. The RPU is typically worn under clothing (e.g., carried in a pocket or purse), but may also be worn in plain sight (e.g., on a belt) if desired. The RPU 16 may be connected (via wired or wireless means 18) to optional secondary wireless link circuitry 19 that allows wireless communication between the RPU and other sources of information (e.g., the general subscriber telephone network) via a secondary wireless link. Note that the optional secondary wireless link circuitry 19 may or may not be contained within the RPU case 70 of FIG. 7.</p> <p data-bbox="583 1125 919 1153">Anderson '783 at 4:30–33</p> <p data-bbox="632 1162 1955 1230">Although the system described here is monaural, it will be clear to those skilled in the art that a binaural (stereo) system can be constructed based on the same principles as the monaural system.</p> <p data-bbox="583 1287 898 1315">Anderson '783 at 5:3–21</p> <p data-bbox="632 1325 1990 1466">During normal operation, speech signals from a nearby talker (and other signals in the ambient audio environment) are picked up by the earpiece wireless microphone 20 and transmitted to the RPU FM receiver 24. The RPU FM receiver 24 output level may be adjusted using the RPU FM receiver 24 volume control. The resulting electrical waveform representing signals from the nearby talker travels through the</p>

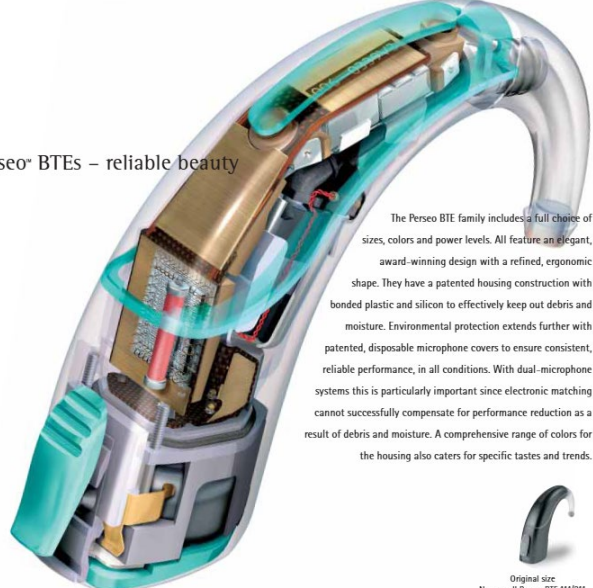
Claims	Prior Art
	<p>DPDT switch 25, which is set to the lower position as shown in FIG. 2, to the signal enhancer 26. The signal enhancer 26 may be, for example, a voice changer device that varies the pitch of a received speech signal according to the setting of pushbutton controls located on the RPU signal enhancer 26. The speech signal's pitch can then be raised or lowered as desired to help compensate for a user's hearing loss relative to a particular talker's voice characteristics. The modified speech signal travels from the RPU signal enhancer 26 to the RPU FM transmitter 27, and finally to the earpiece headset 21 receiver where the signal is converted to acoustic waves heard by the user.</p> <p>Anderson '783 at 6:8–25</p> <p>The user's voice (and other signals in the ambient audio environment) is picked up by the earpiece wireless microphone 20, transmitted to the RPU 23 FM receiver 24, sent through the DPDT switch 25 to the microphone input of the handset 28 (the handset's microphone transducer may be removed) and transmitted on the secondary wireless link 29 to, e.g., the general subscriber telephone network. Signals from the secondary wireless link 29 are received by the handset 28, sent from the speaker output of the handset 28 (the handset's speaker may be removed) via the DPDT switch 25 to the signal enhancer 26, FM transmitter 27 and finally to the earpiece headset 21 receiver where the signals are converted to acoustic waves heard by the user. Note that the RPU 23 may also include straight-forward connections to a "personal digital assistant" computer, or voice-operated devices connected to the FM receiver 24 output and not shown in FIG. 2.</p> <p>Anderson '783 at 21:50–62</p> <p>In the simplest case, the ASIC 845 reads data from the A/D 880, performs a multiplication of the data by a constant (read from the default parameter storage 870) to provide amplification, compares the result with maximum allowable values (read from the default parameter storage 870) and writes the result to the D/A 850. The ASIC 845 may perform other processing (e.g., adaptive volume control) as well. Although none of the remote processor and communication features are available in stand-alone mode, which uses pre-selected hearing aid parameter settings from the default parameter storage 870 in the earpiece, the earpiece provides a useful amplified audio path from the microphone to the speaker.</p> <p>Anderson '783 at 22:11–61</p> <p>Stand-alone earpiece operation: If an earpiece cannot communicate with an RPU for any reason, stand-alone earpiece operation as a CIC hearing aid is generally desirable. Examples include occasions where the RPU must be turned off by pushing the RPU power-off key 78 (e.g., situations where operation of wireless transmitters is prohibited such as aircraft takeoffs and classified military briefings), RPU battery</p>

Claims	Prior Art
	<p>or other RPU electronics failure, misplacement of an RPU, and primary wireless link 17 interference. Note that a microphone located in the earpiece 10 is preferred to a body-worn microphone, since only an earpiece microphone 12 can be used to provide the desired default CIC hearing aid function (earpiece microphones also provide superior audio localization capabilities as the microphone moves with the user's head, and provide reduced wind and clothing noise compared with body-worn microphones). Parameters that control the stand-alone earpiece operation are stored in the earpiece CIC default parameter storage 870, which is accessed by the ASIC 845 to implement the desired stand-alone functions as described earlier. CIC default parameters include, but are not limited to, settings for the nominal and maximum allowable earpiece volume. Other information contained in the default parameter storage 870 includes, but is not limited to, the communication address and L/R (left/right) selection bit. To ensure that the CIC default parameter storage information cannot be changed inadvertently, changes are only allowed when the earpiece is not in use as, for example, when the earpiece battery 860 is being charged. A battery charging detector 840, that is sensitive to the high-level magnetic induction fields commonly used for earpiece battery charging, indicates to the ASIC 845 that the default parameter storage 870 values may be changed upon receipt of an appropriate command from the RPU. Specifically, to change information in the earpiece CIC default parameter storage 870, a user first activates the earpiece battery charging system, e.g., by placing the earpiece in a recharger/storage tray well known in the art. The RPU keyboard 946 is then used to enter a command to the DSP 948 that indicates the user's desire to change information in the default parameter storage 870, followed by the new information to be stored. The DSP 948 causes the RPU to transmit an interrogation with a reserved address code to which all earpieces being charged within the communication range of the RPU are receptive. The information that is to ultimately be loaded into the earpiece default parameter storage 870 is then transmitted by the RPU using the auxiliary bit field of many successive interrogations, said interrogations are received by the earpiece, and the desired information is loaded into the CIC default parameter storage 870 under control of the ASIC 845.</p> <p>Anderson '783 at 22:61–23:35</p> <p>Details of normal system operation: During normal operation of the hearing aid with wireless remote processor, sounds from the ambient environment are picked up by a microphone 12 in the earpiece 10 and transmitted with other information (e.g., data representing the earpiece battery 860 voltage as obtained by an A/D 865 converter and transmitted under control of the earpiece ASIC 845 in the auxiliary bit of many successive replies via the reply data register 825) over a primary two-way wireless link 17 to the RPU 16, where the audio signals are enhanced according to the user's needs by the RPU DSP 948. Signal enhancement is achieved via any of a number of techniques well known in the art that use the capabilities of the general-purpose RPU DSP 948. The enhanced audio signals may be combined with other information (e.g., a synthesized voice, generated within the DSP 948, warning that the aforementioned</p>

Claims	Prior Art
	<p>earpiece battery voltage is low) and transmitted from the RPU 16 over the primary wireless link 17 to the earpiece 10, where they are converted by a speaker 15 to sounds that can only be heard by the user 11. Art auxiliary speaker can be supplied in the RPU 16 as a peripheral 950 device that allows people other than the user to hear the same sounds as those produced by the earpiece speaker 15, if desired (as, e.g., when a person other than the user wishes to use the system's optional cellular telephone capabilities). A different, optional secondary two-way wireless link (shown as part of the RPU in FIG. 9 as secondary wireless link circuitry 944, or connected to the RPU by a wired or other, e.g. infrared, wireless link 18 and shown as optional secondary wireless link circuitry 19 in FIG. 1) can be used for communication between the RPU and a cellular telephone system or other sources of information. An RPU keyboard suitable for covert data entry 72, or voice recognition capabilities in the RPU (implemented, e.g., using any of a number of techniques well known in the art that use the general-purpose RPU DSP 948), can be used to control hearing aid parameters (such as amplification level set by the RPU DSP 948) and telephone dialing functions (e.g., automatic dialing of a predefined telephone number using information stored in the RPU DSP 948) in a manner that can be made imperceptible to a casual observer if desired.</p> <p>Anderson '783 at 23:36–42</p> <p>Binaural, directional noise-canceling hearing aid with hearing protection: Two earpieces (one in each ear) and an RPU can be used in a binaural, directional noise-canceling hearing aid system with hearing protection features. Such a system provides some measure of protection for the residual hearing capability of a hearing-impaired user simultaneous with the hearing aid function.</p> <p>Anderson '783 at 25:29–30</p> <p>The verbal commands are received as audio signals by the earpiece 10 and transmitted to the RPU 16.</p> <p>Anderson '783 at 26:54–65</p> <p>Secondary link audio considerations: In some applications, a user may not wish to transmit signals from the ambient audio environment over the secondary wireless link. In such cases, an auxiliary earpiece microphone which favors the "bone conduction" path over the normal "air path" between the user's mouth and ear can be used. To save space, the auxiliary microphone may be housed in the same miniature case as the primary earpiece microphone. Alternatively, when two earpieces are used, unwanted ambient sources can be "nulled out" in the RPU DSP 948 by digital signal processing techniques well known in the art (e.g., the equivalent of speakerphone circuitry).</p>

Claims	Prior Art
	<p>Anderson '783 at 26:65–67 Note that such binaural signal processing techniques can also provide directivity and background noise cancellation for users</p> <p>Anderson '783 at 27:1–3 having total hearing loss in one ear and partial loss in the other ear, as discussed earlier in connection with binaural hearing aid system operation.</p> <p>A POSITA would have been motivated to combine the teachings of the DHA system with Anderson to meet the claim limitations with a reasonable expectation of success for the reasons described above with respect to claims [1a] and 11.</p> <p>To the extent that the DHA system, alone or in view of Anderson and Nemirovski, does not disclose any elements of this claim limitation, it would have been obvious to a person having ordinary skill in the art to combine, with the teachings of DHA systems, those elements, as known to them or as disclosed by other prior art identified in the cover pleading.</p>
<p>[64b] wherein the first microphone is configured to pick up at least audio input from at least the user,</p>	<p>DHA systems, alone or in view of Anderson, expressly or inherently discloses “wherein the first microphone is configured to pick up at least audio input from at least the user.”</p> <p><i>See</i> above at claim 6 and [50a].</p> <p>A POSITA would have been motivated to combine the teachings of the DHA system with Anderson to meet the claim limitations with a reasonable expectation of success for the reasons described above with respect to claims [1a] and 11.</p> <p>To the extent that the DHA system, alone or in view of Anderson, does not disclose any elements of this claim limitation, it would have been obvious to a person having ordinary skill in the art to combine, with the teachings of DHA systems, those elements, as known to them or as disclosed by other prior art identified in the cover pleading.</p>
<p>[64c] wherein the second microphone is configured to pick up at least environmental audio input associated with an environment in the vicinity of the</p>	<p>DHA systems, alone or in view of Anderson, expressly or inherently discloses “wherein the second microphone is configured to pick up at least environmental audio input associated with an environment in the vicinity of the head-worn wearable electronic device.”</p> <p><i>See</i> above at claim 6 and [50b].</p>

Claims	Prior Art
<p>head-worn wearable electronic device, and</p>	<p>A POSITA would have been motivated to combine the teachings of the DHA system with Anderson to meet the claim limitations with a reasonable expectation of success for the reasons described above with respect to claims [1a] and 11.</p> <p>To the extent that the DHA system, alone or in view of Anderson, does not disclose any elements of this claim limitation, it would have been obvious to a person having ordinary skill in the art to combine, with the teachings of DHA systems, those elements, as known to them or as disclosed by other prior art identified in the cover pleading.</p>
<p>[64d] wherein the processing unit is configured to process the audio input from at least the user via the first microphone to remove at least a portion of the environmental audio input associated with the environment picked up via the second microphone.</p>	<p>DHA systems, alone or in view of Anderson and Nemirovski, expressly or inherently discloses “wherein the processing unit is configured to process the audio input from at least the user via the first microphone to remove at least a portion of the environmental audio input associated with the environment picked up via the second microphone.”</p> <p><i>See above at claim 11.</i></p> <p>A POSITA would have been motivated to combine the teachings of the DHA system with Anderson and Nemirovski to meet the claim limitations with a reasonable expectation of success for the reasons described above with respect to claims [1a] and 11.</p> <p>To the extent that the DHA system, alone or in view of Anderson and Nemirovski, does not disclose any elements of this claim limitation, it would have been obvious to a person having ordinary skill in the art to combine, with the teachings of DHA systems, those elements, as known to them or as disclosed by other prior art identified in the cover pleading.</p>
<p>[65a] A personal audio system as recited in claim 64, wherein the audio input from at least the user picked-up via the first microphone includes at least a voice input, and</p>	<p>DHA systems, alone or in view of Anderson and Nemirovski, expressly or inherently discloses “[a] personal audio system as recited in claim 64, wherein the audio input from at least the user picked-up via the first microphone includes at least a voice input.”</p>

Claims	Prior Art
	<p data-bbox="583 175 1062 207">Phonak – Perseo Brochure at Page 23</p> <div data-bbox="625 224 1260 824">  <p data-bbox="636 383 905 407">Perseo® BTEs – reliable beauty</p> <p data-bbox="968 443 1255 695">The Perseo BTE family includes a full choice of sizes, colors and power levels. All feature an elegant, award-winning design with a refined, ergonomic shape. They have a patented housing construction with bonded plastic and silicon to effectively keep out debris and moisture. Environmental protection extends further with patented, disposable microphone covers to ensure consistent, reliable performance, in all conditions. With dual-microphone systems this is particularly important since electronic matching cannot successfully compensate for performance reduction as a result of debris and moisture. A comprehensive range of colors for the housing also caters for specific tastes and trends.</p> <p data-bbox="1140 800 1247 824">Original size New small Perseo BTE 111/211</p> </div> <p data-bbox="583 846 1625 878">Valente - Effect of type of noise and loudspeaker array on the performance at 399</p> <p data-bbox="636 898 1022 1157">Recently, Phonak introduced its digital hearing aid in which one feature is an adaptive directional microphone. This adaptive microphone reportedly ensures identical sensitivity and phase characteristics of the dual microphones and a mechanism to <i>automatically</i> change the polar pattern from omnidirectional to any directional pattern while compensating for the changes in low-frequency response associated with each pattern. This hearing aid can be programmed</p>

Claims

Prior Art

Phonak – The Hearing Review Article re Perseo (Sep. 5, 2003) at Page 2

Perseo is a 20-channel digital hearing instrument. The default base-program is optimized for speech understanding in quiet environments. A second base-program with different frequency response settings and optional adaptive directional microphones provides additional support in noisy situations. The acoustic environment is continuously monitored and the appropriate base-program is automatically selected. This dual-base automatic function allows for independent fine tuning in one listening situation without compromising performance in the other. The automatic performance is designed to be easily personalized for the individual needs and preferences of the user.

Phonak – Collections 2002 Product List at Page 7



Adaptive digital AudioZoom (dAZ)



Fine-scale Noise Canceller (FNC)

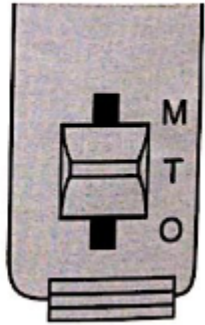
Active and adaptive separation of speech from noise. The Speech-to-Noise Ratio (SNR) can only be optimized if the noise is spatially separated from the speech signal with directional microphone technology. As a pioneer in multi-microphone technologies, Phonak introduced AudioZoom, a feature which significantly improves listening and understanding in noisy situations, enhancing user satisfaction. Now, Phonak goes a step further in the refinement of multi-microphone technology with Adaptive digital AudioZoom: it seeks, finds, and suppresses noise automatically.

Fully automatic, fully adaptive. With its unique spatial processing system, Phonak's Adaptive digital AudioZoom can detect noise not only from directly behind, but anywhere in the near plane. Whether the hearing-impaired person is moving or the location of the noise source changes, the directional characteristic is constantly adapted to the listening situation. In this manner, the Adaptive digital AudioZoom automatically establishes the best possible SNR.

Twenty-band noise control. Claro constantly monitors the Speech-to-Noise Ratio (SNR) in all 20 bands. If a poor SNR is measured in certain bands, the Fine-scale Noise Canceller will reduce gain only in those bands. The result of this precision noise control is a selective and intelligent improvement of the SNR with only a marginal influence on the speech signal.

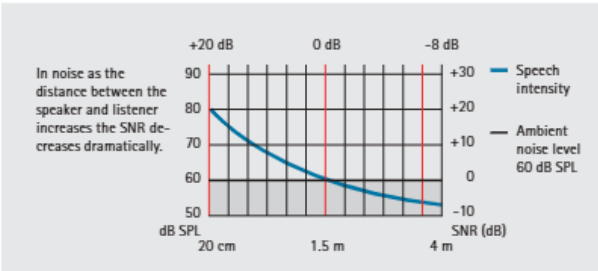
Precise and selective. All frequency bands are separately analysed and controlled. Gain is reduced selectively only in those bands where noise is present. This high-resolution analysis and gain control makes it possible to suppress noise with extremely high accuracy. Neighboring frequency bands without noise remain unaffected.

Intelligent. The magnitude of gain reduction is not constant. In frequency bands with higher SNRs and in bands which are important for speech intelligibility, gain is reduced only slightly. This intelligent control algorithm reduces the masking and distractive effects of background noise. Important speech information is preserved.

Claims	Prior Art
	<p data-bbox="583 175 1444 207">Oticon – DigiFocus Compact User Instructional Brochure at Page 4</p> <div data-bbox="636 228 856 261">Function switch</div> <div data-bbox="636 272 1262 342"> <p>M: Hearing instrument switched on and microphone connected. Normal use.</p> </div> <div data-bbox="636 380 1308 597"> <p>T: Hearing instrument switched on and telecoil connected. Where a telemagnetic loop system has been installed or when you talk on the phone (in many countries), you may then be able to use your hearing instruments with minimum disturbance from background noise.</p> </div> <div data-bbox="636 634 1115 667"> <p>O: Hearing instrument switched off.</p> </div> <div data-bbox="1367 321 1570 643">  </div> <p data-bbox="583 699 932 732">Anderson '783 at 26:54–65</p> <p data-bbox="632 737 1986 987">Secondary link audio considerations: In some applications, a user may not wish to transmit signals from the ambient audio environment over the secondary wireless link. In such cases, an auxiliary earpiece microphone which favors the "bone conduction" path over the normal "air path" between the user's mouth and ear can be used. To save space, the auxiliary microphone may be housed in the same miniature case as the primary earpiece microphone. Alternatively, when two earpieces are used, unwanted ambient sources can be "nulled out" in the RPU DSP 948 by digital signal processing techniques well known in the art (e.g., the equivalent of speakerphone circuitry).</p> <p data-bbox="583 1045 940 1078">Nemirovski '379 at 4:27–41</p> <p data-bbox="632 1083 1990 1403">The communications system 10 includes an ear microphone assembly, or earset 18, having a microphone 20 and a speaker 22 supported by a housing 30 (FIG. 3). The physical arrangement of the earset 18 will be described in more detail below. The microphone 20 is used to detect sounds emanating or coming from the ear (sometimes referred to as air pressure changes within the ear) of the user that result from speech of the user. The microphone 20 converts those detections into an electrical signal. The electrical signal is optionally amplified by an amplifier 24. Next, the amplified electrical signal is received by the communication device 14 via an appropriate interface adapter 26. The interface adapter 26 can be, for example, a jack or other electrical connector received by a corresponding receptacle of the communication device 14.</p>

Claims	Prior Art
	<p>A POSITA would have been motivated to combine the teachings of the DHA system with Anderson and Nemirovski to meet the claim limitations with a reasonable expectation of success for the reasons described above with respect to claims [1a] and 11.</p> <p>To the extent that the DHA system, alone or in view of Anderson and Nemirovski, does not disclose any elements of this claim limitation, it would have been obvious to a person having ordinary skill in the art to combine, with the teachings of DHA systems, those elements, as known to them or as disclosed by other prior art identified in the cover pleading.</p>
<p>[65b] wherein the environmental audio input associated with the environment includes at least ambient noise.</p>	<p>DHA systems, alone or in view of Anderson and Nemirovski, expressly or inherently discloses “wherein the environmental audio input associated with the environment includes at least ambient noise.”</p> <p>Phonak – Perseo Brochure at Page 7</p> <div data-bbox="636 621 1818 1141"> <div>  <div> <p>Convenience and Comfort</p> <p>Perseo is easier to wear and easier to fit</p> <ul style="list-style-type: none"> ■ For true “wear and hear” convenience and comfort in all situations, Perseo automatic operation is managed by Phonak-exclusive PersonalLogic. ■ This breakthrough intelligence for automatic hearing instruments responds to acoustics, the listening situation and personal listening preferences. ■ Perseo is easily programmed for individual needs using the Listening Situation Manager. </div> </div> <div>  <div> <p>Clarity</p> <p>Enhanced performance in multiple listening situations</p> <ul style="list-style-type: none"> ■ The Clarity Component Enhancement system realized in Digital Perception Processing² (DPP²) increases the clarity of speech and environmental sounds. ■ Rapid-tracking Adaptive digital AudioZoom (dAZ) is automatically activated to improve speech understanding in noisy situations. ■ For extra comfort and ease of listening, the Fine-scale Noise Canceler selectively reduces background noise in 20 channels. ■ Optional MicroLink FM provides the ultimate Speech-to-Noise Ratio enhancement. </div> </div> </div>

Claims	Prior Art
	<p data-bbox="583 175 1062 207">Phonak – Perseo Brochure at Page 13</p> <div data-bbox="627 266 1482 308">  <h2 data-bbox="638 266 1482 308">Enhanced sound quality and Speech-to-Noise Ratio</h2> </div> <p data-bbox="726 362 1062 383">Clarity Component Enhancement – DPP²</p> <p data-bbox="726 401 1719 617">Speech clarity depends on both its acoustic properties and the listener's auditory perception. Conventional spectral enhancement systems consider only the acoustics of sound, not how it will be perceived by the listener. Perseo goes well beyond this with DPP², the new generation of Digital Perception Processing. DPP² with Clarity Component Enhancement is individualized for each wearer. It analyzes the acoustic and the psychoacoustic characteristics of incoming sound in 20 critical bands to identify accurately the components that contribute to clarity. It then selectively amplifies them to optimize clarity for the wearer.</p> <p data-bbox="583 638 1024 670">Phonak – Perseo Brochure at p. 14</p> <p data-bbox="638 688 840 716">Rapid-tracking dAZ</p> <p data-bbox="638 742 1056 1122">Phonak pioneered the development of directional microphone systems. Not surprisingly, Perseo dAZ with new soft switching brings even better performance. Perseo dAZ is automatically activated in appropriate situations and immediately reduces noise from beside and behind the wearer. The wearer experiences an immediate reduction of background noise and enhanced speech clarity. Rapid tracking ensures that even fast moving noise sources are controlled.</p>

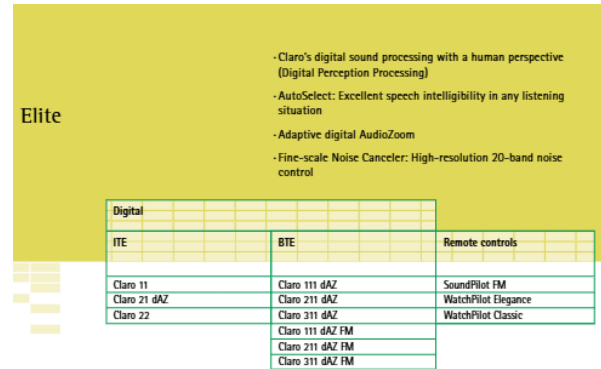
Claims	Prior Art
	<p data-bbox="583 175 1062 207">Phonak – Perseo Brochure at Page 15</p> <p data-bbox="636 228 1358 264">The perfectly integrated communication system</p> <p data-bbox="720 269 831 289">FM MicroLink</p> <p data-bbox="720 313 999 630"> FM technology has allowed thousands of children to participate in standard education settings and to be actively involved in social activities. Phonak pioneered this development and the integration of FM into hearing instruments to produce communication systems that are effective in even the most demanding situations. Since the innovation of miniaturized Phonak MicroLink systems, FM has become an important solution also for adults with a high communication need in work and leisure environments. </p> <p data-bbox="720 667 999 764"> All Perseo BTEs have been designed to incorporate MicroLink with its fully integrated appearance and operation. This integrated system is unique to Phonak. </p> <div data-bbox="1062 313 1656 824"> <p data-bbox="1062 313 1608 332">Speech-to-Noise Ratio (SNR) as a function of distance between speaker and listener</p>  <p data-bbox="1062 634 1432 654">Speech-to-Noise Ratio without FM is distance dependent</p> <div data-bbox="1062 667 1656 824"> <div data-bbox="1062 667 1272 824"> <p data-bbox="1140 675 1199 695">1.5 m</p>  <p data-bbox="1087 805 1243 824">Without FM SNR = 0 dB</p> </div> <div data-bbox="1339 667 1656 824"> <p data-bbox="1482 675 1520 695">4 m</p>  <p data-bbox="1411 805 1587 824">Without FM SNR = -8 dB</p> </div> </div> </div>

Claims	Prior Art
	<p data-bbox="583 175 1062 207">Phonak – Perseo Brochure at Page 23</p> <div data-bbox="625 224 1260 824">  <p data-bbox="634 381 903 406">Perseo® BTEs – reliable beauty</p> <p data-bbox="968 440 1255 691">The Perseo BTE family includes a full choice of sizes, colors and power levels. All feature an elegant, award-winning design with a refined, ergonomic shape. They have a patented housing construction with bonded plastic and silicon to effectively keep out debris and moisture. Environmental protection extends further with patented, disposable microphone covers to ensure consistent, reliable performance, in all conditions. With dual-microphone systems this is particularly important since electronic matching cannot successfully compensate for performance reduction as a result of debris and moisture. A comprehensive range of colors for the housing also caters for specific tastes and trends.</p> <p data-bbox="1140 800 1247 821">Original size New small Perseo BTE 111/211</p> </div>

Claims

Prior Art

Phonak – Collections 2002 Product List at Page 5



Elite Line features:

- Clear, natural sound with precise, fully automatic volume control
- High level of user convenience with automatic hearing program selection
- FM receiver fully integrated in the BTE case (Claro d2z FM)
- Significantly improved speech intelligibility even in difficult listening situations
- Innovative high-tech solutions for Claro system control
- Design innovations for all system components
- State-of-the-art hearing computers for the optimised solution to a wide range of hearing problems

Claims

Prior Art

Phonak – Collections 2002 Product List at Page 7



Adaptive digital AudioZoom (dAZ)



Fine-scale Noise Canceler (FNC)

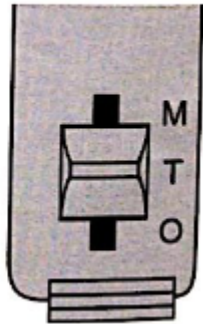
Active and adaptive separation of speech from noise
The Speech-to-Noise Ratio (SNR) can only be optimised if the noise is spatially separated from the speech signal with directional microphone technology. As a pioneer in multi-microphone technologies, Phonak introduced AudioZoom, a feature which significantly improves listening and understanding in noisy situations, enhancing user satisfaction. Now, Phonak goes a step further in the refinement of multi-microphone technology with Adaptive digital AudioZoom: it seeks, finds, and suppresses noise automatically.

Fully automatic, fully adaptive
With its unique spatial processing system, Phonak's Adaptive digital AudioZoom can detect noise not only from directly behind, but anywhere in the near plane. Whether the hearing-impaired person is moving or the location of the noise source changes, the directional characteristic is constantly adapted to the listening situation. In this manner, the Adaptive digital AudioZoom automatically establishes the best possible SNR.

Twenty-band noise control
Claro constantly monitors the Speech-to-Noise Ratio (SNR) in all 20 bands. If a poor SNR is measured in certain bands, the Fine-scale Noise Canceler will reduce gain only in those bands. The result of this precision noise control is a selective and intelligent improvement of the SNR with only a marginal influence on the speech signal.

Precise and selective
All frequency bands are separately analysed and controlled. Gain is reduced selectively only in those bands where noise is present. This high-resolution analysis and gain control makes it possible to suppress noise with extremely high accuracy. Neighboring frequency bands without noise remain unaffected.

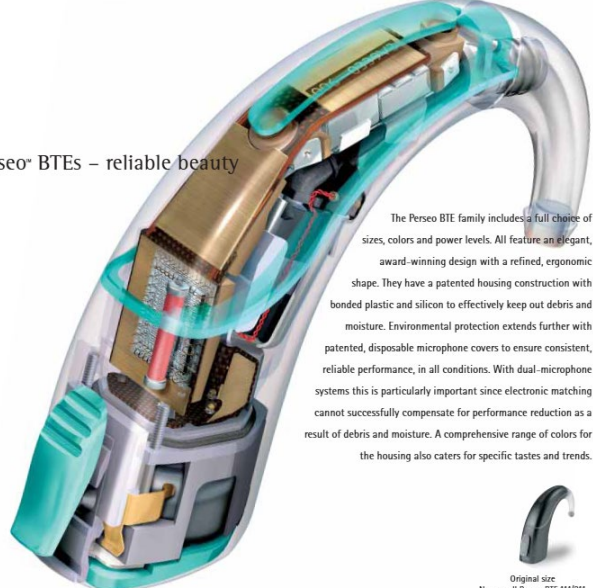
Intelligent
The magnitude of gain reduction is not constant. In frequency bands with higher SNRs and in bands which are important for speech intelligibility, gain is reduced only slightly. This intelligent control algorithm reduces the masking and distractive effects of background noise. Important speech information is preserved.

Claims	Prior Art
	<p data-bbox="583 175 1608 212">Lewis – Speech perception in noise: directional microphones (2004) at Page 431</p> <p data-bbox="636 228 1016 818"> All subjects were fit with digital Phonak Claro 311 dAZ BTE hearing aids bilaterally. All earmolds had select-a-vent (SAV) venting and #13 or 3 mm horn tubing. In addition to the hearing aids, subjects were fit with Phonak Microlink ML8 FM receivers bilaterally. These FM receivers attach to the bottom of a BTE hearing aid and may be utilized in either the “FM only” mode, which attenuates the hearing aid microphone by 20 dB, or in the FM plus hearing aid mode, which allows for FM input and input of environmental sounds via the hearing aid microphones simultaneously at the same output level. The Phonak TX3 HandyMic FM transmitter served as the FM transmitter. This transmitter has three microphone options: (1) Wide Angle, which picks up sounds arriving from all directions around the transmitter microphone equally; (2) Zoom, which provides reduced amplification to signals arriving from the rear (cardioid); and (3) SuperZoom, which provides reduced amplification to signals originating from the rear and the sides (hypercardioid). The </p> <p data-bbox="583 842 1446 880">Oticon – DigiFocus Compact User Instructional Brochure at Page 4</p> <div data-bbox="636 896 858 928"> <p>Function switch</p> </div> <div data-bbox="636 937 1310 1334"> <p>M: Hearing instrument switched on and microphone connected. Normal use.</p> <p>T: Hearing instrument switched on and telecoil connected. Where a telemagnetic loop system has been installed or when you talk on the phone (in many countries), you may then be able to use your hearing instruments with minimum disturbance from background noise.</p> <p>O: Hearing instrument switched off.</p> </div>  <p data-bbox="583 1364 934 1395">Anderson '783 at 26:54–65</p> <p data-bbox="632 1401 1957 1471"> Secondary link audio considerations: In some applications, a user may not wish to transmit signals from the ambient audio environment over the secondary wireless link. In such cases, an auxiliary earpiece </p>

Claims	Prior Art
	<p>microphone which favors the "bone conduction" path over the normal "air path" between the user's mouth and ear can be used. To save space, the auxiliary microphone may be housed in the same miniature case as the primary earpiece microphone. Alternatively, when two earpieces are used, unwanted ambient sources can be "nulled out" in the RPU DSP 948 by digital signal processing techniques well known in the art (e.g., the equivalent of speakerphone circuitry).</p> <p>A POSITA would have been motivated to combine the teachings of the DHA system with Anderson and Nemirovski to meet the claim limitations with a reasonable expectation of success for the reasons described above with respect to claims [1a] and 11.</p> <p>To the extent that the DHA system, alone or in view of Anderson and Nemirovski, does not disclose any elements of this claim limitation, it would have been obvious to a person having ordinary skill in the art to combine, with the teachings of DHA systems, those elements, as known to them or as disclosed by other prior art identified in the cover pleading.</p>
<p>66. A personal audio system as recited in claim 64, wherein the processing unit is configured to process the audio input from at least the user picked-up via the first microphone to remove at least a portion of the environmental audio input picked-up via the second microphone.</p>	<p>DHA systems, alone or in view of Anderson and Nemirovski, expressly or inherently discloses “[a] personal audio system as recited in claim 64, wherein the processing unit is configured to process the audio input from at least the user picked-up via the first microphone to remove at least a portion of the environmental audio input picked-up via the second microphone.”</p> <p><i>See above at claim 11.</i></p> <p>A POSITA would have been motivated to combine the teachings of the DHA system with Anderson and Nemirovski to meet the claim limitations with a reasonable expectation of success for the reasons described above with respect to claims [1a] and 11.</p> <p>To the extent that the DHA system, alone or in view of Anderson and Nemirovski, does not disclose any elements of this claim limitation, it would have been obvious to a person having ordinary skill in the art to combine, with the teachings of DHA systems, those elements, as known to them or as disclosed by other prior art identified in the cover pleading.</p>
<p>[67a] A personal audio system as recited in claim 64, wherein the first microphone is at least partially in the head-worn structure and is positioned so as to be</p>	<p>DHA systems, alone or in view of Anderson and Nemirovski, expressly or inherently discloses “[a] personal audio system as recited in claim 64, wherein the first microphone is at least partially in the head-worn structure and is positioned so as to be directed towards the user of the head-worn structure.”</p> <p><i>See above at claim [50a].</i></p>

Claims	Prior Art
<p>directed towards the user of the head-worn structure, and</p>	<p>A POSITA would have been motivated to combine the teachings of the DHA system with Anderson and Nemirovski to meet the claim limitations with a reasonable expectation of success for the reasons described above with respect to claims [1a] and 11.</p> <p>To the extent that the DHA system, alone or in view of Anderson and Nemirovski, does not disclose any elements of this claim limitation, it would have been obvious to a person having ordinary skill in the art to combine, with the teachings of DHA systems, those elements, as known to them or as disclosed by other prior art identified in the cover pleading.</p>
<p>[67b] wherein the second microphone is at least partially in the head-worn structure so as to be directed away from the user of the head-worn structure.</p>	<p>DHA systems, alone or in view of Anderson and Nemirovski, expressly or inherently discloses “wherein the second microphone is at least partially in the head-worn structure so as to be directed away from the user of the head-worn structure.”</p> <p><i>See above at claim [50b].</i></p> <p>A POSITA would have been motivated to combine the teachings of the DHA system with Anderson and Nemirovski to meet the claim limitations with a reasonable expectation of success for the reasons described above with respect to claims [1a] and 11.</p> <p>To the extent that the DHA system, alone or in view of Anderson and Nemirovski, does not disclose any elements of this claim limitation, it would have been obvious to a person having ordinary skill in the art to combine, with the teachings of DHA systems, those elements, as known to them or as disclosed by other prior art identified in the cover pleading.</p>
<p>68. A personal audio system as recited in claim 67, wherein the processing unit is configured to process the audio input picked-up via at least the first microphone to remove at least a portion of the environmental audio input picked-up via the second microphone.</p>	<p>DHA systems, alone or in view of Anderson and Nemirovski, expressly or inherently discloses “[a] personal audio system as recited in claim 67, wherein the processing unit is configured to process the audio input picked-up via at least the first microphone to remove at least a portion of the environmental audio input picked-up via the second microphone.”</p> <p><i>See above at claim 11.</i></p> <p>A POSITA would have been motivated to combine the teachings of the DHA system with Anderson and Nemirovski to meet the claim limitations with a reasonable expectation of success for the reasons described above with respect to claims [1a] and 11.</p> <p>To the extent that the DHA system, alone or in view of Anderson and Nemirovski, does not disclose any elements of this claim limitation, it would have been obvious to a person having ordinary skill in the art to combine, with the teachings of DHA systems, those elements, as known to them or as disclosed by other prior art identified in the cover pleading.</p>

Claims	Prior Art
<p>[69a] A personal audio system as recited in claim 68, wherein the audio input from at least the user includes at least a voice input, and</p>	<p>DHA systems, alone or in view of Anderson and Nemirovski, expressly or inherently discloses “[a] personal audio system as recited in claim 68, wherein the audio input from at least the user includes at least a voice input.”</p> <p><i>See above at claim [65a].</i></p> <p>A POSITA would have been motivated to combine the teachings of the DHA system with Anderson and Nemirovski to meet the claim limitations with a reasonable expectation of success for the reasons described above with respect to claims [1a] and 11.</p> <p>To the extent that the DHA system, alone or in view of Anderson and Nemirovski, does not disclose any elements of this claim limitation, it would have been obvious to a person having ordinary skill in the art to combine, with the teachings of DHA systems, those elements, as known to them or as disclosed by other prior art identified in the cover pleading.</p>
<p>[69b] wherein the environmental audio input includes at least ambient noise.</p>	<p>DHA systems, alone or in view of Anderson, expressly or inherently discloses “wherein the environmental audio input includes at least ambient noise.”</p> <p><i>See above at claim [65b].</i></p> <p>A POSITA would have been motivated to combine the teachings of the DHA system with Anderson and Nemirovski to meet the claim limitations with a reasonable expectation of success for the reasons described above with respect to claims [1a] and 11.</p> <p>To the extent that the DHA system, alone or in view of Anderson and Nemirovski, does not disclose any elements of this claim limitation, it would have been obvious to a person having ordinary skill in the art to combine, with the teachings of DHA systems, those elements, as known to them or as disclosed by other prior art identified in the cover pleading.</p>
<p>[70a] A personal audio system as recited in claim 24, wherein the head-worn wearable electronic device comprises: a head-worn structure comprising at least two microphones and a processing unit,</p>	<p>DHA systems, alone or in view of Anderson and Nemirovski, expressly or inherently discloses “[a] personal audio system as recited in claim 24, wherein the head-worn wearable electronic device comprises: a head-worn structure comprising at least two microphones and a processing unit.”</p>

Claims	Prior Art
	<p data-bbox="583 175 1062 207">Phonak – Perseo Brochure at Page 23</p> <div data-bbox="625 224 1260 824">  <p data-bbox="636 383 905 407">Perseo® BTEs – reliable beauty</p> <p data-bbox="968 443 1255 695">The Perseo BTE family includes a full choice of sizes, colors and power levels. All feature an elegant, award-winning design with a refined, ergonomic shape. They have a patented housing construction with bonded plastic and silicon to effectively keep out debris and moisture. Environmental protection extends further with patented, disposable microphone covers to ensure consistent, reliable performance, in all conditions. With dual-microphone systems this is particularly important since electronic matching cannot successfully compensate for performance reduction as a result of debris and moisture. A comprehensive range of colors for the housing also caters for specific tastes and trends.</p> <p data-bbox="1140 800 1247 824">Original size New small Perseo BTE 111/211</p> </div> <p data-bbox="583 846 1627 878">Valente - Effect of type of noise and loudspeaker array on the performance at 399</p> <p data-bbox="636 898 1022 1157">Recently, Phonak introduced its digital hearing aid in which one feature is an adaptive directional microphone. This adaptive microphone reportedly ensures identical sensitivity and phase characteristics of the dual microphones and a mechanism to <i>automatically</i> change the polar pattern from omnidirectional to any directional pattern while compensating for the changes in low-frequency response associated with each pattern. This hearing aid can be programmed</p>

Claims	Prior Art
	<p data-bbox="583 175 1514 212">Phonak – The Hearing Review Article re Perseo (Sep. 5, 2003) at Page 2</p> <p data-bbox="636 228 1226 521">Perseo is a 20-channel digital hearing instrument. The default base-program is optimized for speech understanding in quiet environments. A second base-program with different frequency response settings and optional adaptive directional microphones provides additional support in noisy situations. The acoustic environment is continuously monitored and the appropriate base-program is automatically selected. This dual-base automatic function allows for independent fine tuning in one listening situation without compromising performance in the other. The automatic performance is designed to be easily personalized for the individual needs and preferences of the user.</p> <p data-bbox="583 545 1176 583">Phonak – Annual Report (2003–2004) at p. 12</p> <p data-bbox="636 599 1087 989">Within the FM business segment, the Group continued to leverage existing key technology into new applications. With SmartLink, Phonak has achieved another industry first by allowing hearing instrument users to access other communications systems. SmartLink combines a remote control, a transmitter, a directional microphone system and a Bluetooth link to mobile phone all in one single device. Phonak also presented EduLink, an innovative FM device to assist the learning development of children with Auditory Processing Disorder (APD). EduLink was perceived very positively within the medical community and recent classroom tests show favourable results.</p>

Claims

Prior Art

Phonak – Collections 2002 Product List at Page 6

Elite

Digital Perception Processing (DPP)

AutoSelect (AS)

New Phonak sound:
unmatched clarity
Digital Perception Processing (DPP) is based on a cochlear sound analysis across 20 critical bands – exactly like the analysis performed by the healthy human ear. But Clari gets a step further than previous overlapping multi-band filter systems. DPP models the active, load-ress-dependent cochlear signal processing in all 20 critical bands.

20 critical bands:
overlapping and dynamic
The cochlear filters are overlapping, but the degree of overlap is dependent on the frequency, level, and bandwidth of the sound signal. This is precisely how Clari's filters work. Digital Perception Processing computes the cochlear excitation patterns and thus determines sound perception. DPP is the first to mimic the healthy human ear in this way.

Automatic, intelligent program selection
Clari listens and automatically activates the optimal DPP algorithm for the acoustic environment: Program 1 with QuietAdapt for quiet situations or Program 2 with NoiseAdapt and noise control for noisy environments.

In quiet situations
QuietAdapt DPP is designed for clear hearing in quiet surroundings. QuietAdapt controls gain in 20 critical bands to correct for loudness and provide a natural sound quality. By preserving the perception patterns of complex sounds, QuietAdapt assures that low-level signal components are not over-

amplified to the point where they would adversely affect speech intelligibility. This allows Clari to deliver only those sound components that would be naturally audible.

In noise
NoiseAdapt DPP is the ideal program for maximizing audibility in noisy environments. The special NoiseAdapt algorithm works with the Fine-scale Noise Canceller – and if available with the Adaptive digital AudioZoom feature – to optimize speech understanding in noise.

Claims

Prior Art

Phonak – Collections 2002 Product List at Page 7



Adaptive digital AudioZoom (dAZ)

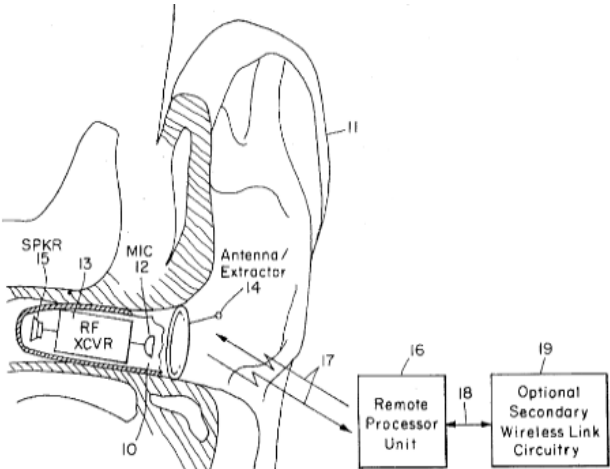


Fine-scale Noise Canceler (FNC)

<p>Active and adaptive separation of speech from noise The Speech-to-Noise Ratio (SNR) can only be optimised if the noise is spatially separated from the speech signal with directional microphone technology. As a pioneer in multi-microphone technologies, Phonak introduced AudioZoom, a feature which significantly improves listening and understanding in noisy situations, enhancing user satisfaction. Now, Phonak goes a step further in the refinement of multi-microphone technology with Adaptive digital AudioZoom: it seeks, finds, and suppresses noise automatically.</p>	<p>Fully automatic, fully adaptive With its unique spatial processing system, Phonak's Adaptive digital AudioZoom can detect noise not only from directly behind, but anywhere in the rear plane. Whether the hearing-impaired person is moving or the location of the noise source changes, the directional characteristic is constantly adapted to the listening situation. In this manner, the Adaptive digital AudioZoom automatically establishes the best possible SNR.</p>	<p>Twenty-band noise control Claro constantly monitors the Speech-to-Noise Ratio (SNR) in all 20 bands. If a poor SNR is measured in certain bands, the Fine-scale Noise Canceler will reduce gain only in those bands. The result of this precision noise control is a selective and intelligent improvement of the SNR with only a marginal influence on the speech signal.</p>	<p>Intelligent The magnitude of gain reduction is not constant. In frequency bands with higher SNRs and in bands which are important for speech intelligibility, gain is reduced only slightly. This intelligent control algorithm reduces the masking and distractive effects of background noise. Important speech information is preserved.</p>
		<p>Precise and selective All frequency bands are separately analysed and controlled. Gain is reduced selectively only in those bands where noise is present. This high-resolution analysis and gain control makes it possible to suppress noise with extremely high accuracy. Neighboring frequency bands without noise remain unaffected.</p>	

Oticon – DigiFocus White Paper (1996) at Page 2

Manufacturers of hearing instruments have worked hard to overcome the troubles of size and battery consumption. With the thickness of digital microchips at only 1/100 of a human hair, this has now been accomplished in the world's first fully digital, at the ear hearing instrument - DigiFocus. Inside DigiFocus there is computing power available which equals that of a PC.

Claims	Prior Art
	<p data-bbox="583 175 894 207">Anderson '783 at FIG. 1</p>  <p data-bbox="995 748 1062 781">FIG. 1</p> <p data-bbox="583 846 919 878">Anderson '783 at 3:13–18</p> <p data-bbox="632 883 1969 1024">FIG. 1 is a block diagram system overview illustrating how an earpiece (comprising a microphone, RF transceiver, dual-use antenna/extractor and speaker) is worn in the ear, communicates via a primary two-way RF link with an RPU, and also communicates via an optional secondary wireless link to a telephone system.</p> <p data-bbox="583 1081 945 1114">Anderson '783 at 3:52–4:14</p> <p data-bbox="632 1118 1986 1443">One preferred embodiment of the invention herein described is shown in block diagram form in FIG. 1. In this figure, the earpiece 10 shown worn in the ear 11 uses a standard completely-in-the-canal (CIC) housing, well known in the art, although many other earpiece housing types known in the art (e.g., behind-the-ear, or BTE) can be used. The earpiece shown in FIG. 1 comprises a microphone 12, RF transceiver 13 with an antenna that doubles as an earpiece extractor 14, and speaker 15. Although the dual-use antenna/extractor arrangement shown in FIG. 1 is often desirable, a separate antenna or antennae can be used in conjunction with many different types of earpiece extractors as appropriate for a particular application. Antennae and extractors can be disguised as ear hair or jewelry (e.g., earrings), or conductive filaments can be permanently implanted in the ear cartilage to act as antennae using a process similar to</p>

Claims	Prior Art
	<p>ear piercing. The earpiece 10 communicates with a remote processor unit (RPU) 16 via a primary two-way RF link 17, although many other wireless link media (e.g., ultrasonic or infrared) may be used instead. The RPU is typically worn under clothing (e.g., carried in a pocket or purse), but may also be worn in plain sight (e.g., on a belt) if desired. The RPU 16 may be connected (via wired or wireless means 18) to optional secondary wireless link circuitry 19 that allows wireless communication between the RPU and other sources of information (e.g., the general subscriber telephone network) via a secondary wireless link. Note that the optional secondary wireless link circuitry 19 may or may not be contained within the RPU case 70 of FIG. 7.</p> <p>Anderson '783 at 5:3–21</p> <p>During normal operation, speech signals from a nearby talker (and other signals in the ambient audio environment) are picked up by the earpiece wireless microphone 20 and transmitted to the RPU FM receiver 24. The RPU FM receiver 24 output level may be adjusted using the RPU FM receiver 24 volume control. The resulting electrical waveform representing signals from the nearby talker travels through the DPDT switch 25, which is set to the lower position as shown in FIG. 2, to the signal enhancer 26. The signal enhancer 26 may be, for example, a voice changer device that varies the pitch of a received speech signal according to the setting of pushbutton controls located on the RPU signal enhancer 26. The speech signal's pitch can then be raised or lowered as desired to help compensate for a user's hearing loss relative to a particular talker's voice characteristics. The modified speech signal travels from the RPU signal enhancer 26 to the RPU FM transmitter 27, and finally to the earpiece headset 21 receiver where the signal is converted to acoustic waves heard by the user.</p> <p>Anderson '783 at 6:8–25</p> <p>The user's voice (and other signals in the ambient audio environment) is picked up by the earpiece wireless microphone 20, transmitted to the RPU 23 FM receiver 24, sent through the DPDT switch 25 to the microphone input of the handset 28 (the handset's microphone transducer may be removed) and transmitted on the secondary wireless link 29 to, e.g., the general subscriber telephone network. Signals from the secondary wireless link 29 are received by the handset 28, sent from the speaker output of the handset 28 (the handset's speaker may be removed) via the DPDT switch 25 to the signal enhancer 26, FM transmitter 27 and finally to the earpiece headset 21 receiver where the signals are converted to acoustic waves heard by the user. Note that the RPU 23 may also include straight-forward connections to a "personal digital assistant" computer, or voice-operated devices connected to the FM receiver 24 output and not shown in FIG. 2.</p>

Claims	Prior Art
	<p data-bbox="583 175 934 207">Anderson '783 at 21:50–62</p> <p data-bbox="632 215 1992 467">In the simplest case, the ASIC 845 reads data from the A/D 880, performs a multiplication of the data by a constant (read from the default parameter storage 870) to provide amplification, compares the result with maximum allowable values (read from the default parameter storage 870) and writes the result to the D/A 850. The ASIC 845 may perform other processing (e.g., adaptive volume control) as well. Although none of the remote processor and communication features are available in stand-alone mode, which uses pre-selected hearing aid parameter settings from the default parameter storage 870 in the earpiece, the earpiece provides a useful amplified audio path from the microphone to the speaker.</p> <p data-bbox="583 524 934 557">Anderson '783 at 22:11–61</p> <p data-bbox="632 565 1992 1474">Stand-alone earpiece operation: If an earpiece cannot communicate with an RPU for any reason, stand-alone earpiece operation as a CIC hearing aid is generally desirable. Examples include occasions where the RPU must be turned off by pushing the RPU power-off key 78 (e.g., situations where operation of wireless transmitters is prohibited such as aircraft takeoffs and classified military briefings), RPU battery or other RPU electronics failure, misplacement of an RPU, and primary wireless link 17 interference. Note that a microphone located in the earpiece 10 is preferred to a body-worn microphone, since only an earpiece microphone 12 can be used to provide the desired default CIC hearing aid function (earpiece microphones also provide superior audio localization capabilities as the microphone moves with the user's head, and provide reduced wind and clothing noise compared with body-worn microphones). Parameters that control the stand-alone earpiece operation are stored in the earpiece CIC default parameter storage 870, which is accessed by the ASIC 845 to implement the desired stand-alone functions as described earlier. CIC default parameters include, but are not limited to, settings for the nominal and maximum allowable earpiece volume. Other information contained in the default parameter storage 870 includes, but is not limited to, the communication address and L/R (left/right) selection bit. To ensure that the CIC default parameter storage information cannot be changed inadvertently, changes are only allowed when the earpiece is not in use as, for example, when the earpiece battery 860 is being charged. A battery charging detector 840, that is sensitive to the high-level magnetic induction fields commonly used for earpiece battery charging, indicates to the ASIC 845 that the default parameter storage 870 values may be changed upon receipt of an appropriate command from the RPU. Specifically, to change information in the earpiece CIC default parameter storage 870, a user first activates the earpiece battery charging system, e.g., by placing the earpiece in a recharger/storage tray well known in the art. The RPU keyboard 946 is then used to enter a command to the DSP 948 that indicates the user's desire to change information in the default parameter storage 870, followed by the new information to be stored. The DSP 948 causes the RPU to transmit an interrogation with a reserved address code to which all earpieces being charged within the communication range of the RPU are receptive. The information that is to ultimately be loaded into</p>

Claims	Prior Art
	<p>the earpiece default parameter storage 870 is then transmitted by the RPU using the auxiliary bit field of many successive interrogations, said interrogations are received by the earpiece, and the desired information is loaded into the CIC default parameter storage 870 under control of the ASIC 845.</p> <p>Anderson '783 at 22:61–23:35</p> <p>Details of normal system operation: During normal operation of the hearing aid with wireless remote processor, sounds from the ambient environment are picked up by a microphone 12 in the earpiece 10 and transmitted with other information (e.g., data representing the earpiece battery 860 voltage as obtained by an A/D 865 converter and transmitted under control of the earpiece ASIC 845 in the auxiliary bit of many successive replies via the reply data register 825) over a primary two-way wireless link 17 to the RPU 16, where the audio signals are enhanced according to the user's needs by the RPU DSP 948. Signal enhancement is achieved via any of a number of techniques well known in the art that use the capabilities of the general-purpose RPU DSP 948. The enhanced audio signals may be combined with other information (e.g., a synthesized voice, generated within the DSP 948, warning that the aforementioned earpiece battery voltage is low) and transmitted from the RPU 16 over the primary wireless link 17 to the earpiece 10, where they are converted by a speaker 15 to sounds that can only be heard by the user 11. Art auxiliary speaker can be supplied in the RPU 16 as a peripheral 950 device that allows people other than the user to hear the same sounds as those produced by the earpiece speaker 15, if desired (as, e.g., when a person other than the user wishes to use the system's optional cellular telephone capabilities). A different, optional secondary two-way wireless link (shown as part of the RPU in FIG. 9 as secondary wireless link circuitry 944, or connected to the RPU by a wired or other, e.g. infrared, wireless link 18 and shown as optional secondary wireless link circuitry 19 in FIG. 1) can be used for communication between the RPU and a cellular telephone system or other sources of information. An RPU keyboard suitable for covert data entry 72, or voice recognition capabilities in the RPU (implemented, e.g., using any of a number of techniques well known in the art that use the general-purpose RPU DSP 948), can be used to control hearing aid parameters (such as amplification level set by the RPU DSP 948) and telephone dialing functions (e.g., automatic dialing of a predefined telephone number using information stored in the RPU DSP 948) in a manner that can be made imperceptible to a casual observer if desired.</p> <p>Anderson '783 at 23:36–42</p> <p>Binaural, directional noise-canceling hearing aid with hearing protection: Two earpieces (one in each ear) and an RPU can be used in a binaural, directional noise-canceling hearing aid system with hearing protection features. Such a system provides some measure of protection for the residual hearing capability</p>

Claims	Prior Art
	<p data-bbox="632 175 1528 212">of a hearing-impaired user simultaneous with the hearing aid function.</p> <p data-bbox="583 266 936 297">Anderson '783 at 25:29–30</p> <p data-bbox="632 303 1929 337">The verbal commands are received as audio signals by the earpiece 10 and transmitted to the RPU 16.</p> <p data-bbox="583 391 936 422">Anderson '783 at 26:54–65</p> <p data-bbox="632 428 1986 683">Secondary link audio considerations: In some applications, a user may not wish to transmit signals from the ambient audio environment over the secondary wireless link. In such cases, an auxiliary earpiece microphone which favors the "bone conduction" path over the normal "air path" between the user's mouth and ear can be used. To save space, the auxiliary microphone may be housed in the same miniature case as the primary earpiece microphone. Alternatively, when two earpieces are used, unwanted ambient sources can be "nulled out" in the RPU DSP 948 by digital signal processing techniques well known in the art (e.g., the equivalent of speakerphone circuitry).</p> <p data-bbox="583 737 936 768">Anderson '783 at 26:65–67</p> <p data-bbox="632 774 1944 841">Note that such binaural signal processing techniques can also provide directivity and background noise cancellation for users</p> <p data-bbox="583 894 905 925">Anderson '783 at 27:1–3</p> <p data-bbox="632 932 1986 1003">having total hearing loss in one ear and partial loss in the other ear, as discussed earlier in connection with binaural hearing aid system operation.</p> <p data-bbox="583 1057 926 1088">Nemirovski '379 at FIG. 5</p>

Claims	Prior Art
	<div data-bbox="583 175 1407 630"> <p style="text-align: center;">FIG. 5</p> </div> <p>Nemirovski '379 at 3:6–23</p> <p>The present invention is directed to a communication system, including a communications device and an ear microphone, or earset, assembly. The communications system allows a user to speak with a remotely located person. The ear microphone assembly, also referred to herein as an earset, includes a microphone and a speaker supported by a housing. The housing is retained by the ear of the user and allows for hands free use of the communication system while carrying on a conversation with the remotely located person. The microphone is directed towards the ear canal of the user and detects sounds emanating or coming out of the ear (or air pressure changes occurring within the ear) to accurately and reliably detect the speech of the user. The position of the microphone may also help to reduce detection of externally generated sounds. The ear microphone assembly allows for separation of a speech input device for a cellular, or wireless, telephone from an RF transceiver of the telephone.</p> <p>Nemirovski '379 at 9:45–60</p> <p>FIG. 5 is a block diagram that illustrates a system 10' for detecting speech of the user. The system 10' includes the earset 18 having the internal microphone 20, the speaker 22 and amplifier 24a as described in greater detail above. In addition, the earset 18 of the system 10' has the external microphone 70. As one skilled in the art will appreciate, the use of the external microphone 70 is optional. The external microphone 70 is disposed on the housing 30 (FIG. 3) and generally points away from the person wearing the earset 18. As a result, the external microphone 70 detects air pressure changes located outside the ear. These air pressure changes are mainly the result of externally generated noise (other people's speech,</p>

Claims	Prior Art
	<p>environmental noise, etc.) and, in part, to the speech of the person wearing the earset 18, when the speech travels from the mouth of the user and through the air to the external microphone 70.</p> <p>A POSITA would have been motivated to combine the teachings of the DHA system with Anderson and Nemirovski to meet the claim limitations with a reasonable expectation of success for the reasons described above with respect to claim [1a] and [51]. For example, like DHA systems, Anderson discloses that its hearing aid system can include an earpiece with two microphones and that certain sounds detected by those microphones can be emphasized over others—for example, “unwanted ambient sources can be ‘nulled out’,” and the system can perform “background noise cancellation.” Anderson, 26:54-67. It would have been obvious for the POSITA to look to Nemirovski to fill in the missing details, as both DHA devices, Anderson and Nemirovski disclose ear-worn devices for hearing impaired individuals. DHA system, Anderson and Nemirovski share various structural and functional similarities that render the inventions readily combinable. For example, both Anderson and Nemirovski disclose ear-worn devices for providing audio output to a user. Both devices include similar components, such as microphones, speakers, wireless transceivers, and control circuitry. Anderson, FIG. 8; Victorian, FIG. 4A. Accordingly, combining the teachings of Anderson and Nemirovski would involve combining known components to yield predictable results.</p> <p>To the extent that the DHA system, alone or in view of Anderson and Nemirovski, does not disclose any elements of this claim limitation, it would have been obvious to a person having ordinary skill in the art to combine, with the teachings of DHA systems, those elements, as known to them or as disclosed by other prior art identified in the cover pleading.</p>
<p>[70b] wherein the at least two microphones are configured to pick up (i) at least audio input from at least the user and (ii) at least environmental audio input associated with environmental sound in the vicinity of the head-worn structure, and</p>	<p>DHA systems, alone or in view of Anderson and Nemirovski, expressly or inherently discloses “wherein the at least two microphones are configured to pick up (i) at least audio input from at least the user and (ii) at least environmental audio input associated with environmental sound in the vicinity of the head-worn structure.”</p>

Claims	Prior Art
	<p>Phonak – Perseo Brochure at Page 7</p> <div data-bbox="636 232 737 332" style="display: inline-block; width: 40px; height: 60px; background-color: #FFD700; margin-right: 10px;"></div> <div data-bbox="758 228 1207 267" style="display: inline-block;">Convenience and Comfort</div> <div data-bbox="758 271 1102 295" style="display: inline-block;">Perseo is easier to wear and easier to fit</div> <div data-bbox="732 354 1224 636" style="display: inline-block; vertical-align: top;"> <ul style="list-style-type: none"> ■ For true "wear and hear" convenience and comfort in all situations, Perseo automatic operation is managed by Phonak-exclusive PersonalLogic. ■ This breakthrough intelligence for automatic hearing instruments responds to acoustics, the listening situation and personal listening preferences. ■ Perseo is easily programmed for individual needs using the Listening Situation Manager. </div> <div data-bbox="1230 232 1331 332" style="display: inline-block; width: 40px; height: 60px; background-color: #00AEEF; margin-right: 10px;"></div> <div data-bbox="1352 228 1470 267" style="display: inline-block;">Clarity</div> <div data-bbox="1352 271 1816 295" style="display: inline-block;">Enhanced performance in multiple listening situations</div> <div data-bbox="1327 354 1820 750" style="display: inline-block; vertical-align: top;"> <ul style="list-style-type: none"> ■ The Clarity Component Enhancement system realized in Digital Perception Processing² (DPP²) increases the clarity of speech and environmental sounds. ■ Rapid-tracking Adaptive digital AudioZoom (dAZ) is automatically activated to improve speech understanding in noisy situations. ■ For extra comfort and ease of listening, the Fine-scale Noise Canceler selectively reduces background noise in 20 channels. ■ Optional MicroLink FM provides the ultimate Speech-to-Noise Ratio enhancement. </div> <p>Phonak – Perseo Brochure at Page 13</p> <div data-bbox="636 862 1480 901" style="text-align: center;"> <h2 style="margin: 0;">Enhanced sound quality and Speech-to-Noise Ratio</h2> </div> <div data-bbox="726 956 1062 977" style="text-align: center;"> <p>Clarity Component Enhancement – DPP²</p> </div> <div data-bbox="726 995 1722 1211" style="text-align: center;"> <p>Speech clarity depends on both its acoustic properties and the listener's auditory perception. Conventional spectral enhancement systems consider only the acoustics of sound, not how it will be perceived by the listener. Perseo goes well beyond this with DPP², the new generation of Digital Perception Processing. DPP² with Clarity Component Enhancement is individualized for each wearer. It analyzes the acoustic and the psychoacoustic characteristics of incoming sound in 20 critical bands to identify accurately the components that contribute to clarity. It then selectively amplifies them to optimize clarity for the wearer.</p> </div> <div data-bbox="1644 820 1827 961" style="text-align: right;">  </div>

Claims

Prior Art

Phonak – Perseo Brochure at p. 14

Rapid-tracking dAZ

Phonak pioneered the development of directional microphone systems. Not surprisingly, Perseo **dAZ** with new soft switching brings even better performance. Perseo **dAZ** is automatically activated in appropriate situations and immediately reduces noise from beside and behind the wearer. The wearer experiences an immediate reduction of background noise and enhanced speech clarity. Rapid tracking ensures that even fast moving noise sources are controlled.

Phonak – Perseo Brochure at Page 15

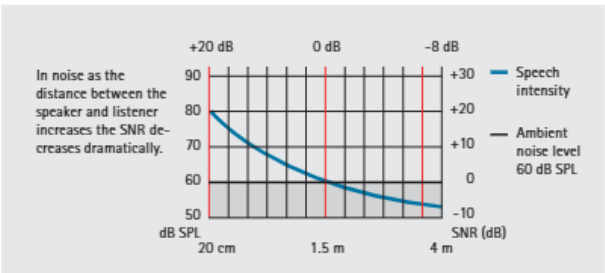
The perfectly integrated communication system

FM MicroLink

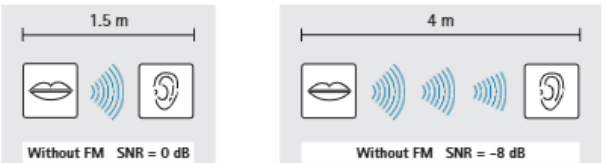
FM technology has allowed thousands of children to participate in standard education settings and to be actively involved in social activities. Phonak pioneered this development and the integration of FM into hearing instruments to produce communication systems that are effective in even the most demanding situations. Since the innovation of miniaturized Phonak **MicroLink** systems, FM has become an important solution also for adults with a high communication need in work and leisure environments.

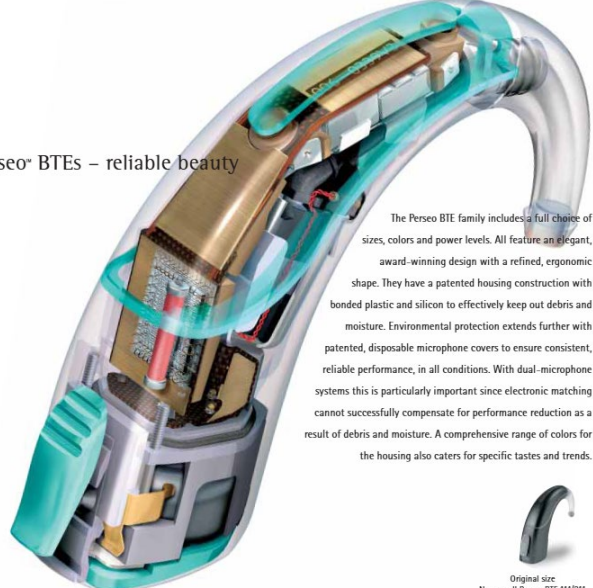
All Perseo BTEs have been designed to incorporate MicroLink with its fully integrated appearance and operation. This integrated system is unique to Phonak.

Speech-to-Noise Ratio (SNR) as a function of distance between speaker and listener



Speech-to-Noise Ratio without FM is distance dependent



Claims	Prior Art
	<p data-bbox="583 175 1062 207">Phonak – Perseo Brochure at Page 23</p> <div data-bbox="625 224 1260 824">  <p data-bbox="636 383 905 407">Perseo® BTEs – reliable beauty</p> <p data-bbox="968 443 1255 695">The Perseo BTE family includes a full choice of sizes, colors and power levels. All feature an elegant, award-winning design with a refined, ergonomic shape. They have a patented housing construction with bonded plastic and silicon to effectively keep out debris and moisture. Environmental protection extends further with patented, disposable microphone covers to ensure consistent, reliable performance, in all conditions. With dual-microphone systems this is particularly important since electronic matching cannot successfully compensate for performance reduction as a result of debris and moisture. A comprehensive range of colors for the housing also caters for specific tastes and trends.</p> <p data-bbox="1140 800 1247 824">Original Size New small Perseo BTE 111/211</p> </div> <p data-bbox="583 846 1627 878">Valente - Effect of type of noise and loudspeaker array on the performance at 399</p> <p data-bbox="636 898 1022 1157">Recently, Phonak introduced its digital hearing aid in which one feature is an adaptive directional microphone. This adaptive microphone reportedly ensures identical sensitivity and phase characteristics of the dual microphones and a mechanism to <i>automatically</i> change the polar pattern from omnidirectional to any directional pattern while compensating for the changes in low-frequency response associated with each pattern. This hearing aid can be programmed</p>

Claims

Prior Art

Phonak – The Hearing Review Article re Perseo (Sep. 5, 2003) at Page 2

Perseo is a 20-channel digital hearing instrument. The default base-program is optimized for speech understanding in quiet environments. A second base-program with different frequency response settings and optional adaptive directional microphones provides additional support in noisy situations. The acoustic environment is continuously monitored and the appropriate base-program is automatically selected. This dual-base automatic function allows for independent fine tuning in one listening situation without compromising performance in the other. The automatic performance is designed to be easily personalized for the individual needs and preferences of the user.

Phonak – Collections 2002 Product List at Page 5

Elite											
- Claro's digital sound processing with a human perspective (Digital Perception Processing)											
- AutoSelect: Excellent speech intelligibility in any listening situation											
- Adaptive digital AudioZoom											
- Fine-scale Noise Canceler: High-resolution 20-band noise control											
Digital											
ITE				BTE				Remote controls			
Claro 11				Claro 111 dAZ				SoundPilot FM			
Claro 21 dAZ				Claro 211 dAZ				WatchPilot Elegance			
Claro 22				Claro 311 dAZ				WatchPilot Classic			
				Claro 111 dAZ FM							
				Claro 211 dAZ FM							
				Claro 311 dAZ FM							

Elite Line features:

- Clear, natural sound with precise, fully automatic volume control
- High level of user convenience with automatic hearing program selection
- FM receiver fully integrated in the BTE case (Claro dAZ FM)
- Significantly improved speech intelligibility even in difficult listening situations
- Innovative high-tech solutions for Claro system control
- State-of-the-art hearing computers for the optimised solution to a wide range of hearing problems
- Design innovations for all system components

Claims

Prior Art

Phonak – Collections 2002 Product List at Page 7



Adaptive digital AudioZoom (dAZ)



Fine-scale Noise Canceler (FNC)

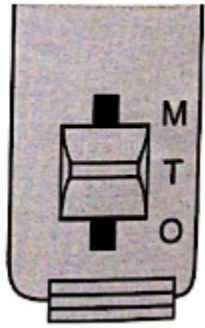
Active and adaptive separation of speech from noise
The Speech-to-Noise Ratio (SNR) can only be optimised if the noise is spatially separated from the speech signal with directional microphone technology. As a pioneer in multi-microphone technologies, Phonak introduced AudioZoom, a feature which significantly improves listening and understanding in noisy situations, enhancing user satisfaction. Now, Phonak goes a step further in the refinement of multi-microphone technology with Adaptive digital AudioZoom: it seeks, finds, and suppresses noise automatically.

Fully automatic, fully adaptive
With its unique spatial processing system, Phonak's Adaptive digital AudioZoom can detect noise not only from directly behind, but anywhere in the near plane. Whether the hearing-impaired person is moving or the location of the noise source changes, the directional characteristic is constantly adapted to the listening situation. In this manner, the Adaptive digital AudioZoom automatically establishes the best possible SNR.

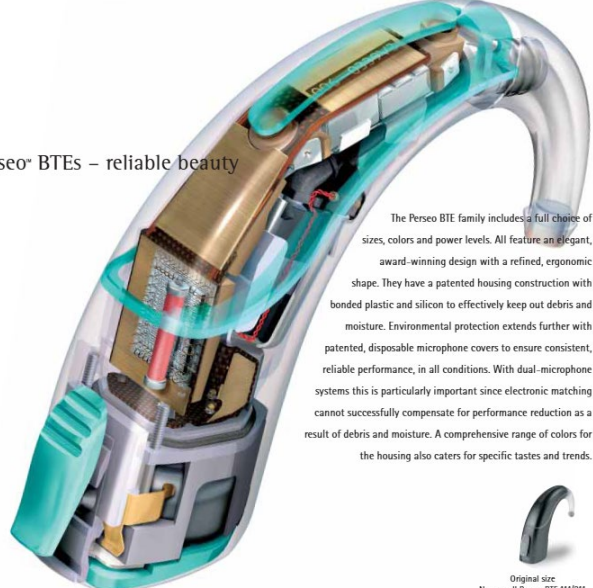
Twenty-band noise control
Claro constantly monitors the Speech-to-Noise Ratio (SNR) in all 20 bands. If a poor SNR is measured in certain bands, the Fine-scale Noise Canceler will reduce gain only in those bands. The result of this precision noise control is a selective and intelligent improvement of the SNR with only a marginal influence on the speech signal.

Precise and selective
All frequency bands are separately analysed and controlled. Gain is reduced selectively only in those bands where noise is present. This high-resolution analysis and gain control makes it possible to suppress noise with extremely high accuracy. Neighboring frequency bands without noise remain unaffected.

Intelligent
The magnitude of gain reduction is not constant. In frequency bands with higher SNRs and in bands which are important for speech intelligibility, gain is reduced only slightly. This intelligent control algorithm reduces the masking and distractive effects of background noise. Important speech information is preserved.

Claims	Prior Art
	<p data-bbox="583 175 1608 212">Lewis – Speech perception in noise: directional microphones (2004) at Page 431</p> <p data-bbox="636 228 1016 818"> All subjects were fit with digital Phonak Claro 311 dAZ BTE hearing aids bilaterally. All earmolds had select-a-vent (SAV) venting and #13 or 3 mm horn tubing. In addition to the hearing aids, subjects were fit with Phonak Microlink ML8 FM receivers bilaterally. These FM receivers attach to the bottom of a BTE hearing aid and may be utilized in either the “FM only” mode, which attenuates the hearing aid microphone by 20 dB, or in the FM plus hearing aid mode, which allows for FM input and input of environmental sounds via the hearing aid microphones simultaneously at the same output level. The Phonak TX3 HandyMic FM transmitter served as the FM transmitter. This transmitter has three microphone options: (1) Wide Angle, which picks up sounds arriving from all directions around the transmitter microphone equally; (2) Zoom, which provides reduced amplification to signals arriving from the rear (cardioid); and (3) SuperZoom, which provides reduced amplification to signals originating from the rear and the sides (hypercardioid). The </p> <p data-bbox="583 842 1446 880">Oticon – DigiFocus Compact User Instructional Brochure at Page 4</p> <div data-bbox="636 894 1310 1333"> <p data-bbox="636 894 858 932">Function switch</p> <p data-bbox="636 938 1262 1008">M: Hearing instrument switched on and microphone connected. Normal use.</p> <p data-bbox="636 1045 1310 1263">T: Hearing instrument switched on and telecoil connected. Where a telemagnetic loop system has been installed or when you talk on the phone (in many countries), you may then be able to use your hearing instruments with minimum disturbance from background noise.</p> <p data-bbox="636 1300 1115 1333">O: Hearing instrument switched off.</p> </div>  <p data-bbox="583 1365 936 1398">Anderson '783 at 26:54–65</p> <p data-bbox="636 1403 1955 1472">Secondary link audio considerations: In some applications, a user may not wish to transmit signals from the ambient audio environment over the secondary wireless link. In such cases, an auxiliary earpiece</p>

Claims	Prior Art
	<p>microphone which favors the "bone conduction" path over the normal "air path" between the user's mouth and ear can be used. To save space, the auxiliary microphone may be housed in the same miniature case as the primary earpiece microphone. Alternatively, when two earpieces are used, unwanted ambient sources can be "nulled out" in the RPU DSP 948 by digital signal processing techniques well known in the art (e.g., the equivalent of speakerphone circuitry).</p> <p>Nemirovski '379 at 4:27–41</p> <p>The communications system 10 includes an ear microphone assembly, or earset 18, having a microphone 20 and a speaker 22 supported by a housing 30 (FIG. 3). The physical arrangement of the earset 18 will be described in more detail below. The microphone 20 is used to detect sounds emanating or coming from the ear (sometimes referred to as air pressure changes within the ear) of the user that result from speech of the user. The microphone 20 converts those detections into an electrical signal. The electrical signal is optionally amplified by an amplifier 24. Next, the amplified electrical signal is received by the communication device 14 via an appropriate interface adapter 26. The interface adapter 26 can be, for example, a jack or other electrical connector received by a corresponding receptacle of the communication device 14.</p> <p>To the extent that the DHA system, alone or in view of Anderson, does not disclose any elements of this claim limitation, it would have been obvious to a person having ordinary skill in the art to combine, with the teachings of DHA systems, those elements, as known to them or as disclosed by other prior art identified in the cover pleading.</p>
<p>[70c] wherein the processing unit is configured to: determine an environmental sound level based on at least the environmental audio input; and</p>	<p>DHA systems, alone or in view of Victorian, expressly or inherently discloses “wherein the processing unit is configured to: determine an environmental sound level based on at least the environmental audio input.”</p>

Claims	Prior Art
	<p data-bbox="583 175 1062 207">Phonak – Perseo Brochure at Page 23</p> <div data-bbox="625 224 1260 824">  <p data-bbox="636 383 905 407">Perseo® BTEs – reliable beauty</p> <p data-bbox="968 443 1255 695">The Perseo BTE family includes a full choice of sizes, colors and power levels. All feature an elegant, award-winning design with a refined, ergonomic shape. They have a patented housing construction with bonded plastic and silicon to effectively keep out debris and moisture. Environmental protection extends further with patented, disposable microphone covers to ensure consistent, reliable performance, in all conditions. With dual-microphone systems this is particularly important since electronic matching cannot successfully compensate for performance reduction as a result of debris and moisture. A comprehensive range of colors for the housing also caters for specific tastes and trends.</p> <p data-bbox="1140 800 1247 824">Original size New small Perseo BTE 111/211</p> </div> <p data-bbox="583 846 1627 878">Valente - Effect of type of noise and loudspeaker array on the performance at 399</p> <p data-bbox="636 898 1020 1157">Recently, Phonak introduced its digital hearing aid in which one feature is an adaptive directional microphone. This adaptive microphone reportedly ensures identical sensitivity and phase characteristics of the dual microphones and a mechanism to <i>automatically</i> change the polar pattern from omnidirectional to any directional pattern while compensating for the changes in low-frequency response associated with each pattern. This hearing aid can be programmed</p>

Claims

Prior Art

Phonak – Collections 2002 Product List at Page 6

Elite



Digital Perception Processing (DPP)



AutoSelect (AS)

New Phonak sound:
unmatched clarity
Digital Perception Processing (DPP) is based on a cochlear sound analysis across 20 critical bands – exactly like the analysis performed by the healthy human ear. But Clari gets a step further than previous overlapping multi-band filter systems. DPP models the active, load-ress-dependent cochlear signal processing in all 20 critical bands.

20 critical bands:
overlapping and dynamic
The cochlear filters are overlapping, but the degree of overlap is dependent on the frequency, level, and bandwidth of the sound signal. This is precisely how Clari's filters work. Digital Perception Processing computes the cochlear excitation patterns and thus determines sound perception. DPP is the first to mimic the healthy human ear in this way.

Automatic, intelligent program selection
Clari listens and automatically activates the optimal DPP algorithm for the acoustic environment: Program 1 with QuietAdapt for quiet situations or Program 2 with NoiseAdapt and noise control for noisy environments.

In quiet situations
QuietAdapt DPP is designed for clear hearing in quiet surroundings. QuietAdapt controls gain in 20 critical bands to correct for loudness and provide a natural sound quality. By preserving the perception patterns of complex sounds, QuietAdapt assures that low-level signal components are not over-

amplified to the point where they would adversely affect speech intelligibility. This allows Clari to deliver only those sound components that would be naturally audible.

In noise
NoiseAdapt DPP is the ideal program for maximizing audibility in noisy environments. The special NoiseAdapt algorithm works with the Fine-scale Noise Canceller – and if available with the Adaptive digital AudioZoom feature – to optimize speech understanding in noise.



Claims

Prior Art

Phonak – Collections 2002 Product List at Page 7



Adaptive digital AudioZoom (dAZ)



Fine-scale Noise Canceler (FNC)

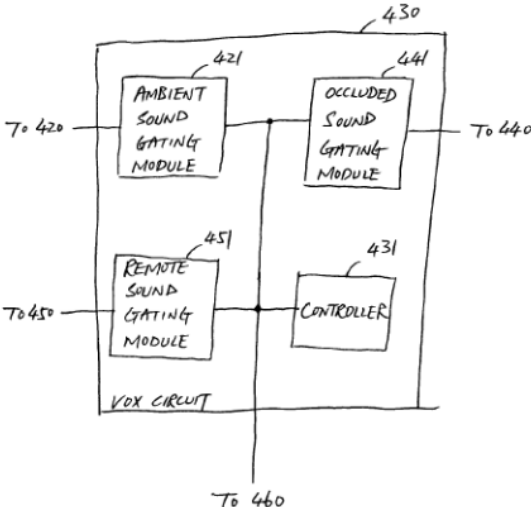
Active and adaptive separation of speech from noise
The Speech-to-Noise Ratio (SNR) can only be optimised if the noise is spatially separated from the speech signal with directional microphone technology. As a pioneer in multi-microphone technologies, Phonak introduced AudioZoom, a feature which significantly improves listening and understanding in noisy situations, enhancing user satisfaction. Now, Phonak goes a step further in the refinement of multi-microphone technology with Adaptive digital AudioZoom: it seeks, finds, and suppresses noise automatically.

Fully automatic, fully adaptive
With its unique spatial processing system, Phonak's Adaptive digital AudioZoom can detect noise not only from directly behind, but anywhere in the near plane. Whether the hearing-impaired person is moving or the location of the noise source changes, the directional characteristic is constantly adapted to the listening situation. In this manner, the Adaptive digital AudioZoom automatically establishes the best possible SNR.

Twenty-band noise control
Claro constantly monitors the Speech-to-Noise Ratio (SNR) in all 20 bands. If a poor SNR is measured in certain bands, the Fine-scale Noise Canceler will reduce gain only in those bands. The result of this precision noise control is a selective and intelligent improvement of the SNR with only a marginal influence on the speech signal.

Intelligent
The magnitude of gain reduction is not constant. In frequency bands with higher SNRs and in bands which are important for speech intelligibility, gain is reduced only slightly. This intelligent control algorithm reduces the masking and distractive effects of background noise. Important speech information is preserved.

Precise and selective
All frequency bands are separately analysed and controlled. Gain is reduced selectively only in those bands where noise is present. This high-resolution analysis and gain control makes it possible to suppress noise with extremely high accuracy. Neighboring frequency bands without noise remain unaffected.

Claims	Prior Art
	<p data-bbox="583 175 915 207">Victorian '313 at FIG. 4B</p>  <p data-bbox="1073 789 1163 821">FIG. 4B</p> <p data-bbox="583 846 894 878">Victorian '313 at [0021]</p> <p data-bbox="632 883 1940 948">FIG. 4B is a block diagram showing one embodiment of portions of a voice operated exchange (VOX) circuit included in the single ear-level device of FIG. 1.</p> <p data-bbox="583 1008 894 1040">Victorian '313 at [0044]</p> <p data-bbox="632 1045 1986 1443">Ear-level device 110 includes an external microphone 420, an ear canal microphone 440, an ear canal speaker 470, a processor 460, a wireless transceiver 450, an antenna 455, a voice operated exchange (VOX) circuit 430, and a battery 480. External microphone 420 picks up ambient sound such that a hearing impaired person wearing ear-level device 110 is not “isolated” when communicating to remote device 190. Thus, ear-level device enables the hearing impaired person to talk to another person either directly (if the other person is nearby) or through a device such as a cell phone. Ear canal microphone 440 detects sound from the ear canal. The sound includes primarily speech of the person wearing ear-level device 110. Ear canal speaker 470 transmits sound received from remote device 190 and/or external microphone 420 to the same ear canal. In one embodiment, ear canal microphone 440 and ear canal speaker 470 are implemented as one physical device. Processor 460 converts the sound picked up from the ear canal to an electrical signal to be transmitted to remote device 190, and converts the signals</p>

Claims	Prior Art
	<p>received from remote device 190 and/or external microphone 420 to a sound audible to the person wearing ear-level device 110. In one embodiment, processor 460 includes one or more of amplification circuitry, filtering circuitry, acoustic feedback reduction circuitry, noise reduction circuitry, and tone control circuitry, among other circuits performing signal processing functions as known in the art. In one embodiment, processor 460 includes a speech recognition module to enhance the audio signal received by and/or transmitted from ear-level device 110. Wireless transceiver 450 and antenna 455 form a telemetry interface supporting telemetry link 195 between ear-level device 110 and remote device 190. Wireless transceiver 450 includes a wireless transmitter and a wireless receiver. The wireless transmitter receives the electrical signal representing the sound picked up from the ear canal from processor 460 and transmits an outgoing wireless signal representing the same sound to remote device 190. The wireless receiver receives an incoming wireless signal representing a remote sound from remote device 190 and transmits an electrical signal representing the remote sound to processor 460. Telemetry link 195 provides for bi-directional communication allowing simultaneous signal transmission, in both directions, between ear-level device 110 and remote device 190. VOX circuit 430 activates a major portion of ear-level device 110 only when sound is detected by at least one of external microphone 420, ear canal microphone 440, and wireless transceiver 450. In one embodiment, VOX circuit 430 includes a voice-controlled switch connecting between battery 480 and the portions of ear-level device 110 that are voice activated. The switch is driven by the sound detected by at least one of external microphone 420, ear canal microphone 440, and wireless transceiver 450. It is turned on when sound detected by external microphone 420 exceeds a predetermined threshold, when sound detected by ear canal microphone 440 exceeds another predetermined threshold, and/or when wireless transceiver 450 receives the incoming wireless signal representing the remote sound. In one embodiment, as discussed with respect to FIG. 4B below, VOX circuit 430 gates or attenuates one or more sounds detected by external microphone 420, ear canal microphone 440, and wireless transceiver 450, to eliminate or reduce echo and ringing caused by the loop formed by ear canal speaker 470 and ear canal microphone 440. Battery 480 supplies the power needed for the operation of ear-level device 110. In one embodiment, battery 480 is a rechargeable battery that can be recharged without being taken out from ear-level device 110.</p> <p>Victorian '313 at [0045]</p> <p>FIG. 4B is a block diagram showing one embodiment of portions of VOX circuit 430. In this embodiment, VOX circuit 430 includes an ambient sound gating module 421 that blanks or attenuates the sound detected by external microphone 420, an occluded sound gating module 441 that blanks or attenuates the sound detected by ear canal microphone 440, a remote sound gating module 451 that blanks or attenuates the sound detected by wireless transceiver 450, and controller 431 to coordinate the states of the gating modules. Because ear canal microphone 440 is capable of detecting the sound delivered to the</p>

Claims	Prior Art
	<p>ear canal from ear canal speaker 470, there is a need to avoid echo and/or ringing. For example, a remote sound from remote device 190 delivered to the ear canal may be detected as an occluded sound and transmitted back to remote device 190, causing echo. The gating modules function in a coordinated manner to eliminate or reduce such echo and/or ringing caused by breaking, or applying resistance into, the audio loop formed by ear canal speaker 470 and ear canal microphone 440.</p> <p>Victorian '313 at [0047] Ambient sound gating module 421 includes a voice-driven switch or attenuator driven by the output of a comparator that compares the amplitude of a detected ambient sound with a predetermined threshold for the ambient sound. When the amplitude of the detected ambient sound exceeds the threshold for the ambient sound, and when controller 431 permits, ambient sound gating module 421 is turned on. When the amplitude of the detected ambient sound drops below the threshold for the ambient sound, ambient sound gating module 421 is turned off. In one embodiment, after the amplitude of the detected ambient sound drops below the threshold for the ambient sound, ambient sound gating module 421 is turned off after a predetermined delay.</p> <p>A POSITA would have been motivated to combine the teachings of the DHA system with Victorian to meet the claim limitations with a reasonable expectation of success for the reasons described above with respect to claims [1a] and [70a].</p> <p>To the extent that the DHA system, alone or in view of Victorian, does not disclose any elements of this claim limitation, it would have been obvious to a person having ordinary skill in the art to combine, with the teachings of DHA systems, those elements, as known to them or as disclosed by other prior art identified in the cover pleading.</p>
<p>[70d] amplify at least a portion of the audio input from the first microphone, with the amplification being at least in part dependent on the environmental sound level.</p>	<p>DHA systems, alone or in view of Victorian, expressly or inherently discloses “amplify at least a portion of the audio input from the first microphone, with the amplification being at least in part dependent on the environmental sound level.”</p> <p><i>See above at claim [51d].</i></p> <p>To the extent that the DHA system, alone or in view of Anderson, does not disclose any elements of this claim limitation, it would have been obvious to a person having ordinary skill in the art to combine, with the teachings of DHA systems, those elements, as known to them or as disclosed by other prior art identified in the cover pleading.</p>

Claims	Prior Art
<p>[71a] A personal audio system as recited in claim 24, wherein the head-worn wearable electronic device comprises at least a first microphone, a second microphone and a processing unit,</p>	<p>DHA systems, alone or in view of Anderson and Nemirovski, expressly or inherently discloses “[a] personal audio system as recited in claim 24, wherein the head-worn wearable electronic device comprises at least a first microphone, a second microphone and a processing unit.”</p> <p><i>See above at claim [70a].</i></p> <p>A POSITA would have been motivated to combine the teachings of the DHA system with Anderson and Nemirovski to meet the claim limitations with a reasonable expectation of success for the reasons described above with respect to claims [1a] and [70a].</p> <p>To the extent that the DHA system, alone or in view of Anderson, and Nemirovski, does not disclose any elements of this claim limitation, it would have been obvious to a person having ordinary skill in the art to combine, with the teachings of DHA systems, those elements, as known to them or as disclosed by other prior art identified in the cover pleading.</p>
<p>[71b] wherein the first microphone is configured to pick up at least audio sound from at least the user, wherein the second microphone is configured to pick up at least ambient noise associated with an environment in the vicinity of the head-worn wearable electronic device, and</p>	<p>DHA systems, alone or in view of Anderson and Nemirovski, expressly or inherently discloses “wherein the first microphone is configured to pick up at least audio sound from at least the user, wherein the second microphone is configured to pick up at least ambient noise associated with an environment in the vicinity of the head-worn wearable electronic device.”</p> <p><i>See above at claim [70b].</i></p> <p>A POSITA would have been motivated to combine the teachings of the DHA system with Anderson and Nemirovski to meet the claim limitations with a reasonable expectation of success for the reasons described above with respect to claims [1a] and [70a].</p> <p>To the extent that the DHA system, alone or in view of Anderson, and Victorian, does not disclose any elements of this claim limitation, it would have been obvious to a person having ordinary skill in the art to combine, with the teachings of DHA systems, those elements, as known to them or as disclosed by other prior art identified in the cover pleading.</p>
<p>[71c] wherein the processing unit is configured to: determine an ambient noise level based on at least the ambient noise; and</p>	<p>DHA systems, alone or in view of Victorian, expressly or inherently discloses “wherein the processing unit is configured to: determine an ambient noise level based on at least the ambient noise.”</p> <p><i>See above at claim [51c].</i></p> <p>A POSITA would have been motivated to combine the teachings of the DHA system with Anderson and Victorian to meet the claim limitations with a reasonable expectation of success for the reasons described above with respect to claims [1a] and [70a].</p>

Claims	Prior Art
	<p>To the extent that the DHA system, alone or in view of Victorian, does not disclose any elements of this claim limitation, it would have been obvious to a person having ordinary skill in the art to combine, with the teachings of DHA systems, those elements, as known to them or as disclosed by other prior art identified in the cover pleading.</p>
<p>[71d] modify at least a portion of the audio sound picked-up via the first microphone, with the modification being at least in part dependent on the ambient noise level.</p>	<p>DHA systems, alone or in view of Victorian, expressly or inherently discloses “modify at least a portion of the audio sound picked-up via the first microphone, with the modification being at least in part dependent on the ambient noise level.”</p> <p><i>See</i> above at claim [51d].</p> <p>A POSITA would have been motivated to combine the teachings of the DHA system with Victorian to meet the claim limitations with a reasonable expectation of success for the reasons described above with respect to claims [1a] and [70a].</p> <p>To the extent that the DHA system, alone or in view of Victorian, does not disclose any elements of this claim limitation, it would have been obvious to a person having ordinary skill in the art to combine, with the teachings of DHA systems, those elements, as known to them or as disclosed by other prior art identified in the cover pleading.</p>
<p>[72a] A personal audio system as recited in claim 24, wherein the head-worn wearable electronic device comprises at least a first microphone, a second microphone and a processing unit,</p>	<p>DHA systems, alone or in view of Anderson and Nemirovski, expressly or inherently discloses “modify at least a portion of the audio sound picked-up via the first microphone, with the modification being at least in part dependent on the ambient noise level.”</p> <p><i>See</i> above at claim [70a].</p> <p>A POSITA would have been motivated to combine the teachings of the DHA system with Anderson Nemirovski to meet the claim limitations with a reasonable expectation of success for the reasons described above with respect to claim [1a], [26a], and [70a].</p> <p>To the extent that the DHA system, alone or in view of Anderson and Nemirovski, does not disclose any elements of this claim limitation, it would have been obvious to a person having ordinary skill in the art to combine, with the teachings of DHA systems, those elements, as known to them or as disclosed by other prior art identified in the cover pleading.</p>
<p>[72b] wherein the first microphone is configured to pick up at least audio sound from at least the user, wherein the second microphone is</p>	<p>DHA systems, alone or in view of Anderson and Nemirovski, expressly or inherently discloses “wherein the first microphone is configured to pick up at least audio sound from at least the user, wherein the second</p>

Claims	Prior Art
<p>configured to pick up at least ambient noise associated with an environment in the vicinity of the head-worn wearable electronic device,</p>	<p>microphone is configured to pick up at least ambient noise associated with an environment in the vicinity of the head-worn wearable electronic device.”</p> <p><i>See above at claim [70b].</i></p> <p>A POSITA would have been motivated to combine the teachings of the DHA system with Anderson and Nemirovski to meet the claim limitations with a reasonable expectation of success for the reasons described above with respect to claim [1a], [26a], and [70a].</p> <p>To the extent that the DHA system, alone or in view of Anderson and Nemirovski, does not disclose any elements of this claim limitation, it would have been obvious to a person having ordinary skill in the art to combine, with the teachings of DHA systems, those elements, as known to them or as disclosed by other prior art identified in the cover pleading.</p>
<p>[72c] wherein the processing unit is configured to: determine an ambient noise level based on at least the ambient noise; and</p>	<p>DHA systems, alone or in view of Anderson and Victorian, expressly or inherently discloses “wherein the processing unit is configured to: determine an ambient noise level based on at least the ambient noise.”</p> <p><i>See above at claim [51c].</i></p> <p>A POSITA would have been motivated to combine the teachings of the DHA system with Anderson, Jannard and Victorian to meet the claim limitations with a reasonable expectation of success for the reasons described above with respect to claim [1a], [26a], and [70a].</p> <p>To the extent that the DHA system, alone or in view of Anderson and Victorian, does not disclose any elements of this claim limitation, it would have been obvious to a person having ordinary skill in the art to combine, with the teachings of DHA systems, those elements, as known to them or as disclosed by other prior art identified in the cover pleading.</p>
<p>[72d] modify at least a portion of the audio sound picked-up via the first microphone, with the modification being at least in part dependent on the ambient noise level,</p>	<p>DHA systems, alone or in view of Anderson and Victorian, expressly or inherently discloses “modify at least a portion of the audio sound picked-up via the first microphone, with the modification being at least in part dependent on the ambient noise level.”</p> <p><i>See above at claim [51d].</i></p> <p>A POSITA would have been motivated to combine the teachings of the DHA system with Anderson, Jannard and Victorian to meet the claim limitations with a reasonable expectation of success for the reasons described above with respect to claim [1a], [26a], and [70a].</p> <p>To the extent that the DHA system, alone or in view of Anderson and Victorian, does not disclose any elements of this claim limitation, it would have been obvious to a person having ordinary skill in the art to</p>

Claims	Prior Art
	<p>combine, with the teachings of DHA systems, those elements, as known to them or as disclosed by other prior art identified in the cover pleading.</p>
<p>[72e] wherein the head-worn wearable electronic device comprises at least a rechargeable battery and at least one conductive element,</p>	<p>DHA systems, alone or in view of Anderson and Jannard, expressly or inherently discloses “wherein the head-worn wearable electronic device comprises at least a rechargeable battery and at least one conductive element.”</p> <p><i>See above at claim [26a].</i></p> <p>A POSITA would have been motivated to combine the teachings of the DHA system with Anderson, Jannard and Victorian to meet the claim limitations with a reasonable expectation of success for the reasons described above with respect to claim [1a], [26a], and [70a].</p> <p>To the extent that the DHA system, alone or in view of Anderson and Victorian, does not disclose any elements of this claim limitation, it would have been obvious to a person having ordinary skill in the art to combine, with the teachings of DHA systems, those elements, as known to them or as disclosed by other prior art identified in the cover pleading.</p>
<p>[72f] wherein the rechargeable battery is able to be electrically coupled via the at least one conductive element to a charging apparatus and charged when the head-worn wearable electronic device is placed within the charging apparatus,</p>	<p>DHA systems, alone or in view of Anderson and Jannard, expressly or inherently discloses “wherein the rechargeable battery is able to be electrically coupled via the at least one conductive element to a charging apparatus and charged when the head-worn wearable electronic device is placed within the charging apparatus.”</p> <p><i>See above at claim [26b].</i></p> <p>A POSITA would have been motivated to combine the teachings of the DHA system with Anderson, Jannard and Victorian to meet the claim limitations with a reasonable expectation of success for the reasons described above with respect to claim [1a], [26a], and [70a].</p> <p>To the extent that the DHA system, alone or in view of Anderson and Victorian, does not disclose any elements of this claim limitation, it would have been obvious to a person having ordinary skill in the art to combine, with the teachings of DHA systems, those elements, as known to them or as disclosed by other prior art identified in the cover pleading.</p>
<p>[72g] wherein the mobile communication device is configured to wirelessly communicate with the head-worn</p>	<p>DHA systems, alone or in view of Anderson, expressly or inherently discloses “wherein the mobile communication device is configured to wirelessly communicate with the head-worn wearable electronic device over a wireless network.”</p> <p><i>See above at claim 27.</i></p>

Claims	Prior Art
wearable electronic device over a wireless network, and	<p>A POSITA would have been motivated to combine the teachings of the DHA system with Anderson, Jannard and Victorian to meet the claim limitations with a reasonable expectation of success for the reasons described above with respect to claim [1a], [26a], and [70a].</p> <p>To the extent that the DHA system, alone or in view of Anderson and Victorian, does not disclose any elements of this claim limitation, it would have been obvious to a person having ordinary skill in the art to combine, with the teachings of DHA systems, those elements, as known to them or as disclosed by other prior art identified in the cover pleading.</p>
[72h] wherein the mobile communication device is or includes a mobile phone,	<p>DHA systems, alone or in view of Anderson, expressly or inherently discloses “wherein the mobile communication device is or includes a mobile phone.”</p> <p><i>See above at claim 14.</i></p> <p>A POSITA would have been motivated to combine the teachings of the DHA system with Anderson, Jannard and Victorian to meet the claim limitations with a reasonable expectation of success for the reasons described above with respect to claim [1a], [26a], and [70a].</p> <p>To the extent that the DHA system, alone or in view of Anderson and Victorian, does not disclose any elements of this claim limitation, it would have been obvious to a person having ordinary skill in the art to combine, with the teachings of DHA systems, those elements, as known to them or as disclosed by other prior art identified in the cover pleading.</p>
[72i] wherein the head-worn wearable electronic device is configured to receive a voice input via the first microphone, determine a requested operation from the voice input using voice recognition, and control an operation by the head-worn wearable electronic device based on the requested operation.	<p>DHA systems, alone or in view of Anderson and Jannard or Seshadri , expressly or inherently discloses “wherein the head-worn wearable electronic device is configured to receive a voice input via the first microphone, determine a requested operation from the voice input using voice recognition, and control an operation by the head-worn wearable electronic device based on the requested operation.”</p> <p><i>See above at claim [53b].</i></p> <p>A POSITA would have been motivated to combine the teachings of the DHA system with Anderson and Jannard or Seshadri to meet the claim limitations with a reasonable expectation of success for the reasons described above with respect to claim [1a], [26a], [53b] and [70a].</p> <p>To the extent that the DHA system, alone or in view of Anderson and Victorian, does not disclose any elements of this claim limitation, it would have been obvious to a person having ordinary skill in the art to combine, with the teachings of DHA systems, those elements, as known to them or as disclosed by other prior art identified in the cover pleading.</p>

Claims	Prior Art
<p>[73a] A personal audio system as recited in claim 72, wherein the head-worn wearable electronic device comprises: an infrared emitter configured to at least direct infrared signals towards the user when the head-worn wearable electronic device is being worn;</p>	<p>DHA systems, alone or in view of Anderson and Abreu or Luo or Irvin, expressly or inherently discloses “[a] personal audio system as recited in claim 72, wherein the head-worn wearable electronic device comprises: an infrared emitter configured to at least direct infrared signals towards the user when the head-worn wearable electronic device is being worn.”</p> <p><i>See above at claim [54a].</i></p> <p>A POSITA would have been motivated to combine the teachings of the DHA system with Anderson, Jannard and Victorian to meet the claim limitations with a reasonable expectation of success for the reasons described above with respect to claim [1a], [26a], and [70a].</p> <p>To the extent that the DHA system, alone or in view of Anderson, Jannard and Victorian, does not disclose any elements of this claim limitation, it would have been obvious to a person having ordinary skill in the art to combine, with the teachings of DHA systems, those elements, as known to them or as disclosed by other prior art identified in the cover pleading.</p>
<p>[73b] an infrared detector configured to detect infrared radiation reflected back from the user when the head-worn wearable electronic device is being worn; and</p>	<p>DHA systems, alone or in view of Anderson and Abreu or Luo or Irvin, expressly or inherently discloses “an infrared detector configured to detect infrared radiation reflected back from the user when the head-worn wearable electronic device is being worn.”</p> <p><i>See above at claim [54b].</i></p> <p>A POSITA would have been motivated to combine the teachings of the DHA system with Anderson and Abreu or Luo or Irvin to meet the claim limitations with a reasonable expectation of success for the reasons described above with respect to claim [1a] and [23a].</p> <p>To the extent that the DHA system, alone or in view of Anderson and Abreu or Luo or Irvin, does not disclose any elements of this claim limitation, it would have been obvious to a person having ordinary skill in the art to combine, with the teachings of DHA systems, those elements, as known to them or as disclosed by other prior art identified in the cover pleading.</p>
<p>[73c] electronic circuitry configured to facilitate determining whether the head-worn wearable electronic device is being worn by the user based on at least the detected infrared radiation.</p>	<p>DHA systems, alone or in view of Anderson and Abreu or Luo or Irvin, expressly or inherently discloses “electronic circuitry configured to facilitate determining whether the head-worn wearable electronic device is being worn by the user based on at least the detected infrared radiation.”</p> <p><i>See above at claim 55.</i></p> <p>A POSITA would have been motivated to combine the teachings of the DHA system with Anderson and Abreu or Luo or Irvin to meet the claim limitations with a reasonable expectation of success for the reasons described above with respect to claim [1a] and [23a].</p>

Claims	Prior Art
	<p>To the extent that the DHA system, alone or in view of Anderson and Abreu or Luo or Irvin, does not disclose any elements of this claim limitation, it would have been obvious to a person having ordinary skill in the art to combine, with the teachings of DHA systems, those elements, as known to them or as disclosed by other prior art identified in the cover pleading.</p>
<p>74. A personal audio system as recited in claim 73, wherein the electrical coupling of the at least one conductive element with the charging apparatus is supported by at least a magnetic element.</p>	<p>DHA systems, alone or in view of Anderson and Jannard, expressly or inherently discloses “[a] personal audio system as recited in claim 73, wherein the electrical coupling of the at least one conductive element with the charging apparatus is supported by at least a magnetic element.”</p> <p><i>See</i> above at claim [26c].</p> <p>A POSITA would have been motivated to combine the teachings of the DHA system with Anderson and Jannard to meet the claim limitations with a reasonable expectation of success for the reasons described above with respect to claim [1a] and [26].</p> <p>To the extent that the DHA system, alone or in view of Anderson and Jannard, does not disclose any elements of this claim limitation, it would have been obvious to a person having ordinary skill in the art to combine, with the teachings of DHA systems, those elements, as known to them or as disclosed by other prior art identified in the cover pleading.</p>
<p>[75pre] A personal audio system operable to at least generate customized audio output, comprising: a pair of speakers, each of the speakers being configured to be wearable proximate to an ear of a user; and</p>	<p>DHA systems, alone or in view of Anderson and/or Enzmann, expressly or inherently discloses “[a] personal audio system operable to at least generate customized audio output, comprising: a pair of speakers, each of the speakers being configured to be wearable proximate to an ear of a user.”</p> <p><i>See</i> above at claim [1pre], [24pre], [1a], and [24a].</p> <p>For example, Anderson discloses a hearing aid with an earpiece in communication with a remote processor unit. <i>See, e.g.,</i> Anderson ’783 at 1:7–10; FIG. 1.</p> <p>In some embodiments, such as shown in FIG. 2, the earpiece belongs to a head-worn headset 21. <i>See</i> 4:15–39; FIG. 2.</p>

Claims	Prior Art
	<div data-bbox="594 180 1533 755" data-label="Diagram"> </div> <p data-bbox="1041 797 1142 829">FIG. 2</p> <p data-bbox="583 867 1948 971">The speaker is wearable proximate to the ear of the user when it is included in the headset 21. This system generates a customized audio output for a user of the hearing aid, for example, by enhancing audio signals according to a user's needs. <i>See, e.g.</i>, Anderson '783 at 1:55–63.</p> <p data-bbox="583 992 1982 1203">Anderson '783 also expressly or inherently discloses “a pair of speakers, each of the speakers being configured to be wearable proximate to an ear of a user.” For example, as described for element [24a], “head-worn” means “eyewear or headwear.” For example, as explained above at claim 24, a POSITA would have understood from Anderson's disclosures that headset 21 can include a pair of speakers, one for each ear, as shown above; with each speaker wearable near (and thus proximate to) its respective ear, <i>i.e.</i>, when the headset is worn.</p> <p data-bbox="583 1224 1974 1333">As described for element [1a], Anderson further discloses that “a binaural (stereo) system can be constructed based on the same principles as the monaural system.” <i>See id.</i> at 4:30-33. Accordingly, the POSITA would have understood Anderson to disclose a pair of speakers (<i>i.e.</i>, one for each of the user's two ears).</p> <p data-bbox="583 1354 1986 1458">Anderson illustrates headset 21 and earpiece 22 within the context of a “limited-feature embodiment,” shown in FIG. 2. <i>See</i> Anderson at 4:15-29. However, a POSITA would have understood headset 21 and earpiece 22 and/or the DHA Systems to be combinable with the “full-featured” embodiment of FIG. 1. <i>See id.</i> at 3:52-58.</p>

Claims	Prior Art
	<p>This is because Anderson characterizes FIG. 1 and FIG. 2 as showing two preferred embodiments of the same invention; explains that the headset earpiece shown in FIG. 2 is “functionally similar” to the in-ear earpiece 10 of FIG. 1; and further explains that “many other earpiece housing types” can be used in place of the CIC earpiece housing shown in FIG. 1. <i>See id.</i> at 3:52-58, 4:15-35. A POSITA would have understood from these disclosures that the FIG. 1 system could be flexibly implemented using various housings for the earpiece, including headset 21.</p> <p>A POSITA would have been motivated to use the headset instead of the CIC housing or DHA System housing to realize advantages in cost, simplicity, and commercial availability described by Anderson. Even using headset 21 in a full-featured application, it would provide certain advantages over the in-ear embodiment. Unlike some CIC and/or DHA System housings, the headset did not need to be customized to fit in a user’s ear canal. Additionally, there would have been comfort benefits: many users do not like supporting objects with their ears, a problem addressed by replacing a CIC housing or DHA System housing with a headset. And using a larger headset rather than a small in-ear device would make it easier and less expensive to include the processing and communication capabilities, and would reduce the likelihood of misplacing the device. Anderson provides specific guidance for adapting the headset using “techniques [that] are well known in the art,” and notes that schematics and service manuals for the headset were readily available, demonstrating that the POSITA would have had a reasonable expectation of success in making the combination. <i>See</i> Anderson at 4:26-5:60.</p> <p>To the extent that the DHA system, alone or in view of Anderson and/or Enzmann, does not disclose any elements of this claim limitation, it would have been obvious to a person having ordinary skill in the art to combine, with the teachings of DHA systems, those elements, as known to them or as disclosed by other prior art identified in the cover pleading.</p>
<p>[75a] a mobile communication device at least for supporting audio output, the mobile communication device configured to store at least a hearing characteristic associated with the user, the stored hearing characteristic associated with the user at least characterizes hearing impairment regarding the user's ability to hear,</p>	<p>DHA systems, alone or in view of Anderson and/or Enzmann, expressly or inherently discloses “a mobile communication device at least for supporting audio output, the mobile communication device configured to store at least a hearing characteristic associated with the user, the stored hearing characteristic associated with the user at least characterizes hearing impairment regarding the user's ability to hear.”</p> <p><i>See</i> above at claim [1b].</p> <p>To the extent that the DHA system, alone or in view of Anderson and/or Enzmann, does not disclose any elements of this claim limitation, it would have been obvious to a person having ordinary skill in the art to</p>

Claims	Prior Art
	<p>combine, with the teachings of DHA systems, those elements, as known to them or as disclosed by other prior art identified in the cover pleading.</p>
<p>[75b] the mobile communication device configured to produce customized audio output signals for the user based on at least the stored hearing characteristics associated with the user such that the customized audio output signals at least compensate for the hearing impairment regarding the user's ability to hear,</p>	<p>DHA systems, alone or in view of Anderson, expressly or inherently discloses “the mobile communication device configured to produce customized audio output signals for the user based on at least the stored hearing characteristics associated with the user such that the customized audio output signals at least compensate for the hearing impairment regarding the user's ability to hear.”</p> <p><i>See above at claim [1c].</i></p> <p>To the extent that the DHA system, alone or in view of Anderson, does not disclose any elements of this claim limitation, it would have been obvious to a person having ordinary skill in the art to combine, with the teachings of DHA systems, those elements, as known to them or as disclosed by other prior art identified in the cover pleading.</p>
<p>[75c] the mobile communication device configured to produce the customized audio output signals using at least signal amplification, and the mobile communication device configured to transmit the customized audio output signals to the pair of speakers,</p>	<p>DHA systems, alone or in view of Anderson, expressly or inherently discloses “the mobile communication device configured to produce the customized audio output signals using at least signal amplification, and the mobile communication device configured to transmit the customized audio output signals to the pair of speakers.”</p> <p><i>See above at claims [1d]; [1e]; [24e].</i></p> <p>To the extent that the DHA system, alone or in view of Anderson, does not disclose any elements of this claim limitation, it would have been obvious to a person having ordinary skill in the art to combine, with the teachings of DHA systems, those elements, as known to them or as disclosed by other prior art identified in the cover pleading.</p>
<p>[75d] wherein each of the speakers is configured to receive at least a portion of the customized audio output signals transmitted by the mobile communication device, and</p>	<p>DHA systems, alone or in view of Anderson, expressly or inherently discloses “wherein each of the speakers is configured to receive at least a portion of the customized audio output signals transmitted by the mobile communication device.”</p> <p><i>See above at claims [1f]; [24f].</i></p> <p>To the extent that the DHA system, alone or in view of Anderson, does not disclose any elements of this claim limitation, it would have been obvious to a person having ordinary skill in the art to combine, with the teachings of DHA systems, those elements, as known to them or as disclosed by other prior art identified in the cover pleading.</p>

Claims	Prior Art
<p>[75e] wherein each of the speakers is configured to generate at least a portion of the customized audio output for the user in accordance with at least the portion of the customized audio output signals being received.</p>	<p>DHA systems, alone or in view of Anderson, expressly or inherently discloses “wherein each of the speakers is configured to generate at least a portion of the customized audio output for the user in accordance with at least the portion of the customized audio output signals being received.”</p> <p><i>See above at claims [1g]; [24g].</i></p> <p>To the extent that the DHA system, alone or in view of Anderson, does not disclose any elements of this claim limitation, it would have been obvious to a person having ordinary skill in the art to combine, with the teachings of DHA systems, those elements, as known to them or as disclosed by other prior art identified in the cover pleading.</p>
<p>76. A personal audio system as recited in claim 75, wherein each of the speakers being configured to be at least partially inserted into an ear of the user.</p>	<p>DHA systems, alone or in view of Anderson, expressly or inherently discloses “[a] personal audio system as recited in claim 75, wherein each of the speakers being configured to be at least partially inserted into an ear of the user.”</p> <p>Phonak – Perseo Brochure at p. 8</p> 

Claims	Prior Art
	<p data-bbox="583 175 1627 207">Valente - Effect of type of noise and loudspeaker array on the performance at 398</p> <p data-bbox="636 228 1211 358">Differences in performance between omnidirectional and directional microphones were evaluated between two loudspeaker conditions (single loudspeaker at 180°; diffuse using eight loudspeakers set 45° apart) and two types of noise (steady-state HINT noise; R-Space™ restaurant noise). Twenty-five participants were fit bilaterally with Phonak Perseo hearing aids using the manufacturer's recommended procedure. After wearing the hearing aids for one week, the</p> <p data-bbox="583 383 945 415">Anderson '783 at 3:52–4:14</p> <p data-bbox="630 420 1988 1036">One preferred embodiment of the invention herein described is shown in block diagram form in FIG. 1. In this figure, the earpiece 10 shown worn in the ear 11 uses a standard completely-in-the-canal (CIC) housing, well known in the art, although many other earpiece housing types known in the art (e.g., behind-the-ear, or BTE) can be used. The earpiece shown in FIG. 1 comprises a microphone 12, RF transceiver 13 with an antenna that doubles as an earpiece extractor 14, and speaker 15. Although the dual-use antenna/extractor arrangement shown in FIG. 1 is often desirable, a separate antenna or antennae can be used in conjunction with many different types of earpiece extractors as appropriate for a particular application. Antennae and extractors can be disguised as ear hair or jewelry (e.g., earrings), or conductive filaments can be permanently implanted in the ear cartilage to act as antennae using a process similar to ear piercing. The earpiece 10 communicates with a remote processor unit (RPU) 16 via a primary two-way RF link 17, although many other wireless link media (e.g., ultrasonic or infrared) may be used instead. The RPU is typically worn under clothing (e.g., carried in a pocket or purse), but may also be worn in plain sight (e.g., on a belt) if desired. The RPU 16 may be connected (via wired or wireless means 18) to optional secondary wireless link circuitry 19 that allows wireless communication between the RPU and other sources of information (e.g., the general subscriber telephone network) via a secondary wireless link. Note that the optional secondary wireless link circuitry 19 may or may not be contained within the RPU case 70 of FIG. 7.</p> <p data-bbox="583 1094 1967 1235">To the extent that the DHA system, alone or in view of Anderson, does not disclose any elements of this claim limitation, it would have been obvious to a person having ordinary skill in the art to combine, with the teachings of DHA systems, those elements, as known to them or as disclosed by other prior art identified in the cover pleading.</p>